

GE Healthcare Class C eService Notes are copyright (c) to General Electric Company Inc with all rights reserved. Class C eService Notes issued by GE Healthcare are proprietary to GE Healthcare Systems. Disclosure to persons not employed by GE Healthcare or not in possession of a Advanced Service Agreement with GE Healthcare is prohibited.
Duplication and distribution of any class of eService Notes, in any form by entities other than GE Healthcare is strictly prohibited.

MI&CT SERVICE NOTE SN06110805a

Class C Service information enclosed.

Date: November 5, 2006

Subject: VCT Gantry Torque labels

Effectivity: All CT Products produced after August 2005.

Details:

CT Products began shipping with torque labels on the following components and systems

- Slipping mounting brackets (Lightpseed 5.x – 7.x)
- 48V power supply bracket (Lightspeed 7.x & Discovery VCT only)R
- Rotating base (Lightspeed 7.x systems & Discovery VCT only).

The torque labels were added to the CT gantry as an enhancement to assist the FE in providing easy and quick access to required hardware information. Early versions of the labels may not remain adhesive to the non-machined casting surfaces and could come loose. These labels could detach and fall into the gantry or slipping causing failures. There are no replacement labels. All information can be found in the appropriate service manual. These labels were previously mentioned in service note SN05120915a. SN05120915a requests removing 5131694 series labels due to incorrect information.

GE Healthcare Class C eService Notes are copyright (c) to General Electric Company Inc with all rights reserved. Class C eService Notes issued by GE Healthcare are proprietary to GE Healthcare Systems. Disclosure to persons not employed by GE Healthcare or not in possession of a Advanced Service Agreement with GE Healthcare is prohibited.
Duplication and distribution of any class of eService Notes, in any form by entities other than GE Healthcare is strictly prohibited.

Resolution:

1) 48V POWER SUPPLY LABELS -- LIGHTSPEED 7.X AND DISCOVERY VCT ONLY

- Remove label with the following part numbers on 48V power supply bracket. Part number is shown in the lower right corner of the label.
 - o 5131694-3
 - o 5155151-3
- The following label part numbers have the correct adhesive and **DO NOT** require removal.
 - o 5177970-3

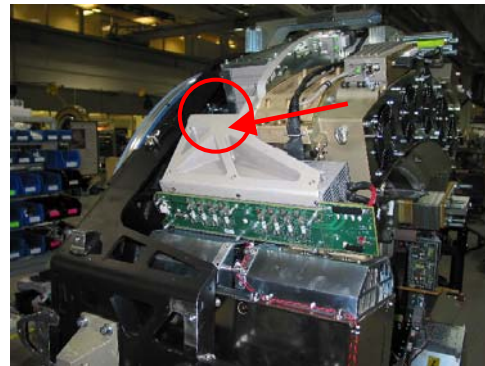


Figure 1: label 5131694-3 or 5155151-3 (48V PWR SUPPLY to DAS BRKT)

2) DETECTOR BOTTOM RAIL -- LIGHTSPEED 7.X & DISCOVERY VCT ONLY

- Remove label with the following part numbers on rotating base near detector collimator. Part number is shown in the lower right corner of the label.
 - o 5131694-6
 - o 5155151-6
- The following label part numbers have the correct adhesive and **DO NOT** require removal.
 - o 5177970-6



Figure 2: label 5131694-6 or 5155151-6 (ROT BASE into DET BOTTOM RAIL)

Duplication and distribution of any class of eService Notes, in any form by entities other than GE Healthcare is strictly prohibited.

3) COLLIMATOR LABEL - LIGHTSPEED 7.X & DISCOVERY VCT ONLY

- Remove label with the following part numbers on the rotating base near the tube collimator. Part number is shown in the lower right corner of the label.
 - o 5131694-7
 - o 5155151-7
- The following label part numbers have the correct adhesive and **DO NOT** require removal.
 - o 5177970-7

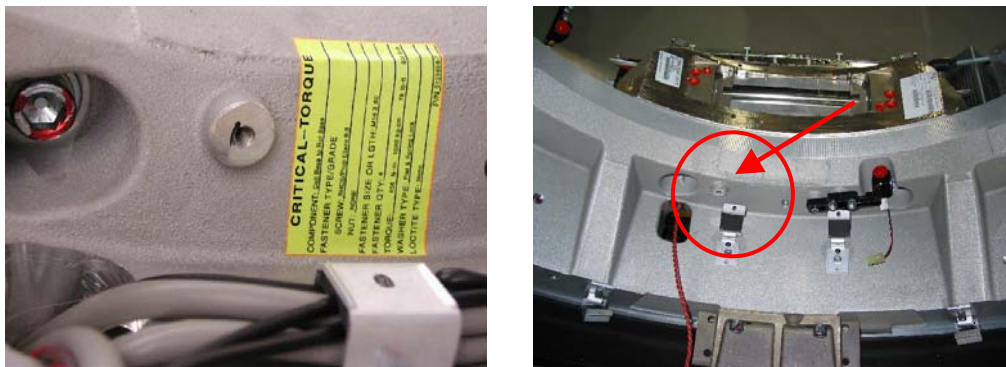
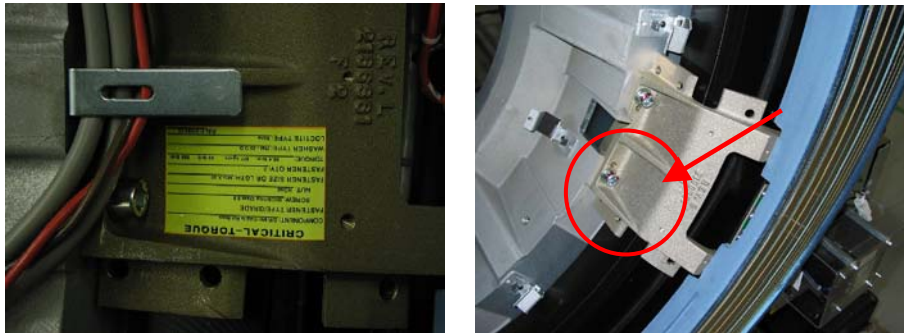


Figure 3: label 5131694-7 or 5155151-7 (COLL BASE to ROT BASE)

4) SLIPRING BRACKET LABELS -- ALL CT and CT/PET SYSTEMS

- Remove slipring brackets labels with the following part numbers. Part number is shown in the lower right corner of the label.
 - o 5131694-14
 - o 5131694-15
 - o 5155151-14
 - o 5155151-15
- The following label part numbers have the correct adhesive and **DO NOT** require removal.
 - o 5177970-14
 - o 5177970-15



**Figure 4: label 5131694-15 or 5155151-15 (S/R MNT CASTING to ROT BASE),
5131694-14 or 5155151-14 not shown**

Shawn Ray
Hardware Quality Team
MI&CT Engineering