

GE Healthcare MI&CT Class C eService Notes are copyright (c) protected to General Electric Company Inc with all rights reserved. Class C eService Notes issued by GE Healthcare are proprietary to GE Healthcare Systems. Disclosure to persons not employed by GE Healthcare or not in possession of a Advanced Service Agreement with GE Healthcare is prohibited.
Duplication and distribution of any class of eService Notes, in any form by entities other than GE Healthcare is strictly prohibited.

MI&CT ESERVICE NOTE SN07021310a

Class C information is Property of GE.

This short-term informational note is invalid 6 months after the issue date.

Date: February 21, 2007

Subject: Global Table Leg Bearings Issue

Effectivity: VCT 64 and 32, RT16 and PET Discovery VCT

Details: It was found that the grease was not correctly applied to the leg bearings of the Global Table manufactured before November 2006. This may cause a serious damage on the leg bearings.

Resolution: Check and grease the bearings according to the following instructions.
This will be also released as FMI later (FMI# 26812).

PRELIMINARY REQUIREMENTS

Consumables:

Qty	Description	Part #
1	Grease Polyex (30cc)	2347076 or
	Shell Alvania Grease #2 (60cc)	2188551
1	Penetrating Spray Lubricant (CRC5-56, WD-40 or equivalent.)	-
1	Kim-Wipe or Waste cloth	-

Tools:

Qty	Description	Part #
1	Grease Gun (for Polyex 2347076)	46-297931P1

PROCEDURES

- 1). Raise the table to maximum height.
- 2). Move the IMS and cradle to OUT limit position.
- 3). Remove the front base cover.
- 4). Remove power from the table circuit breaker.

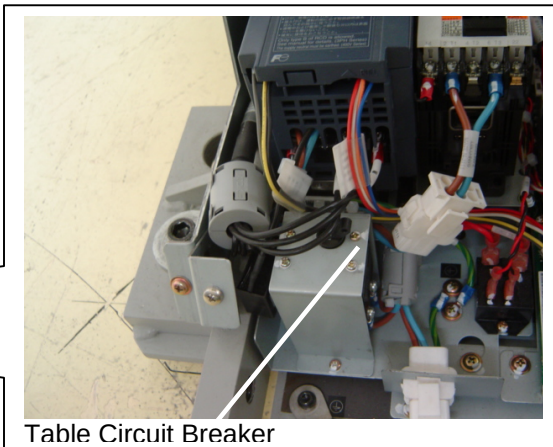
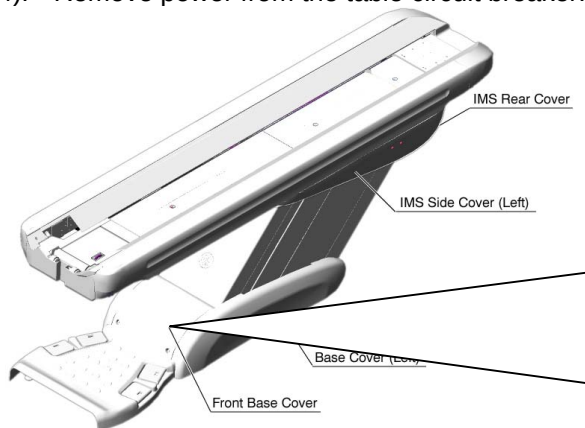


Figure - 1 Table Circuit Breaker

GE Healthcare MI&CT Class C eService Notes are copyright (c) protected to General Electric Company Inc with all rights reserved. Class C eService Notes issued by GE Healthcare are proprietary to GE Healthcare Systems. Disclosure to persons not employed by GE Healthcare or not in possession of a Advanced Service Agreement with GE Healthcare is prohibited.

Duplication and distribution of any class of eService Notes, in any form by entities other than GE Healthcare is strictly prohibited.

5). Remove the following table covers:

- IMS Side Cover (Left/Right)
- IMS Rear Cover
- Base Cover (Left/Right)


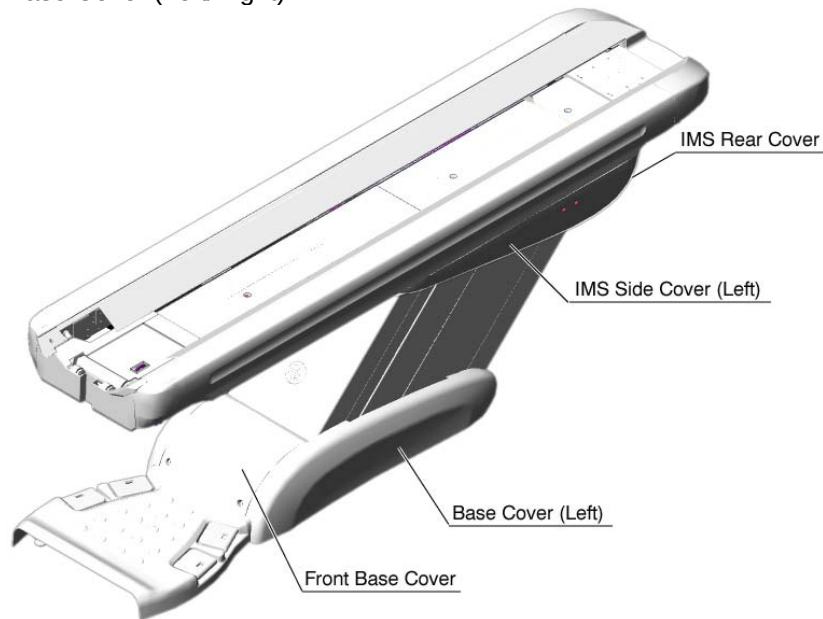


Figure - 2 Table Covers Removal

6). Check if the grease is correctly applied to all eight (8) leg bearings.
If it is OK, go to step 18 to complete this Service Note.
If it is NOT OK, proceed to the next steps.

Grease is applied

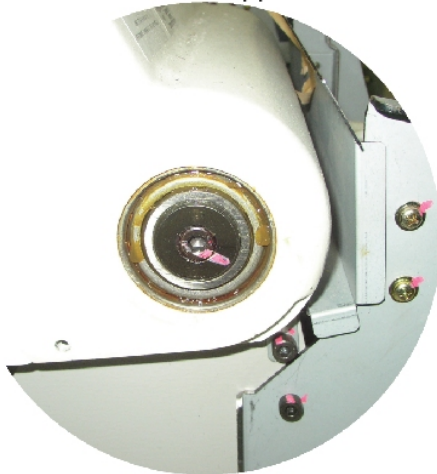

Grease is not applied



Figure - 3 Table Leg Bearing

Duplication and distribution of any class of eService Notes, in any form by entities other than GE Healthcare is strictly prohibited.

7). Remove the following table covers:

- Front Side Plate (Left/Right)
- Rear Side Plate (Left/Right)
- Cover Plate (Left/Right)


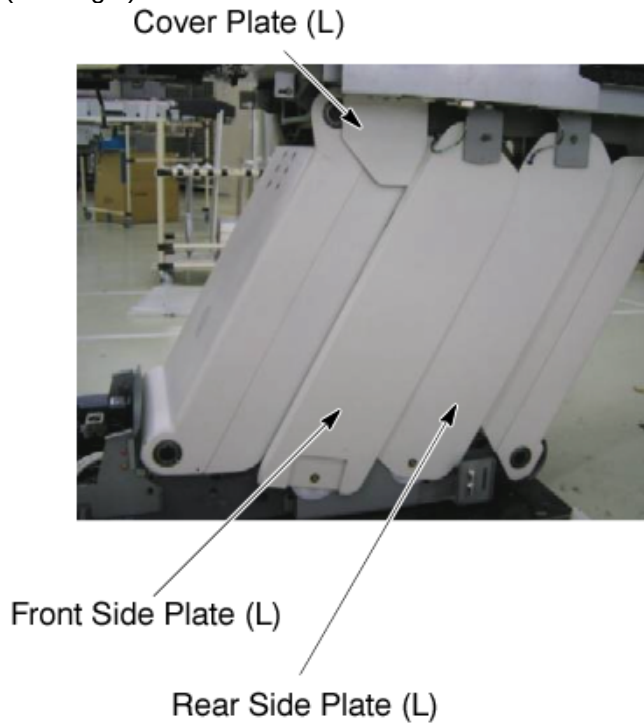


Figure - 4 Side Plates and Cover Plates Removal

8). Cover the Inverter, Power Supply and GTCB to prevent oils contamination.

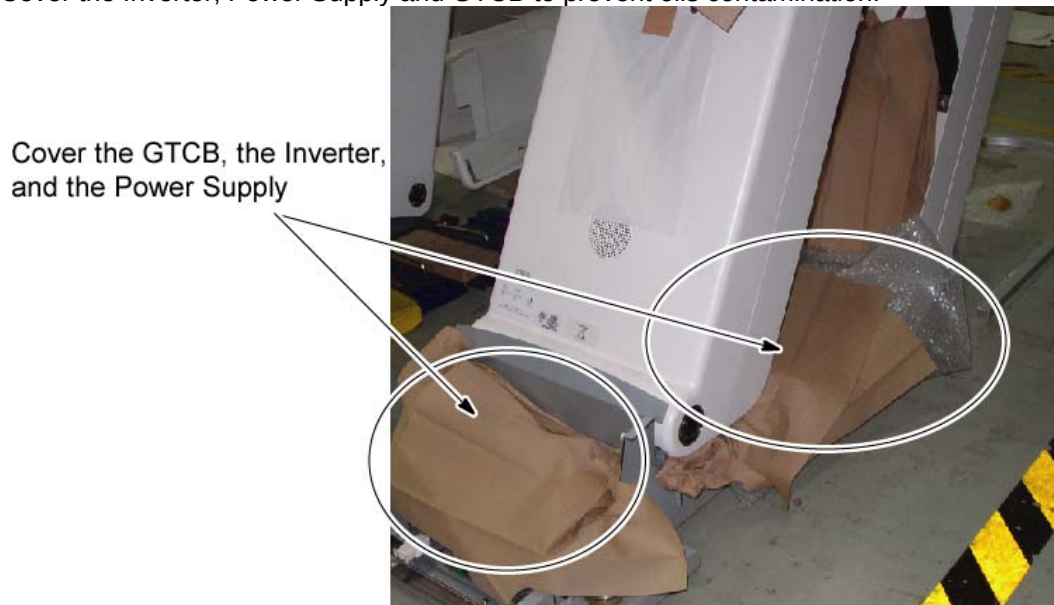

GE Healthcare MI&CT Class C eService Notes are copyright (c) protected to General Electric Company Inc with all rights reserved. Class C eService Notes issued by GE Healthcare are proprietary to GE Healthcare Systems. Disclosure to persons not employed by GE Healthcare or not in possession of a Advanced Service Agreement with GE Healthcare is prohibited.

Duplication and distribution of any class of eService Notes, in any form by entities other than GE Healthcare is strictly prohibited.

Figure - 5 Inverter, Power Supply and GTCB Protection

- 9). Cram Kim-Wipe (or waste cloth) into each joint gap between the table leg and the leg bearing to prevent the oil dripping.

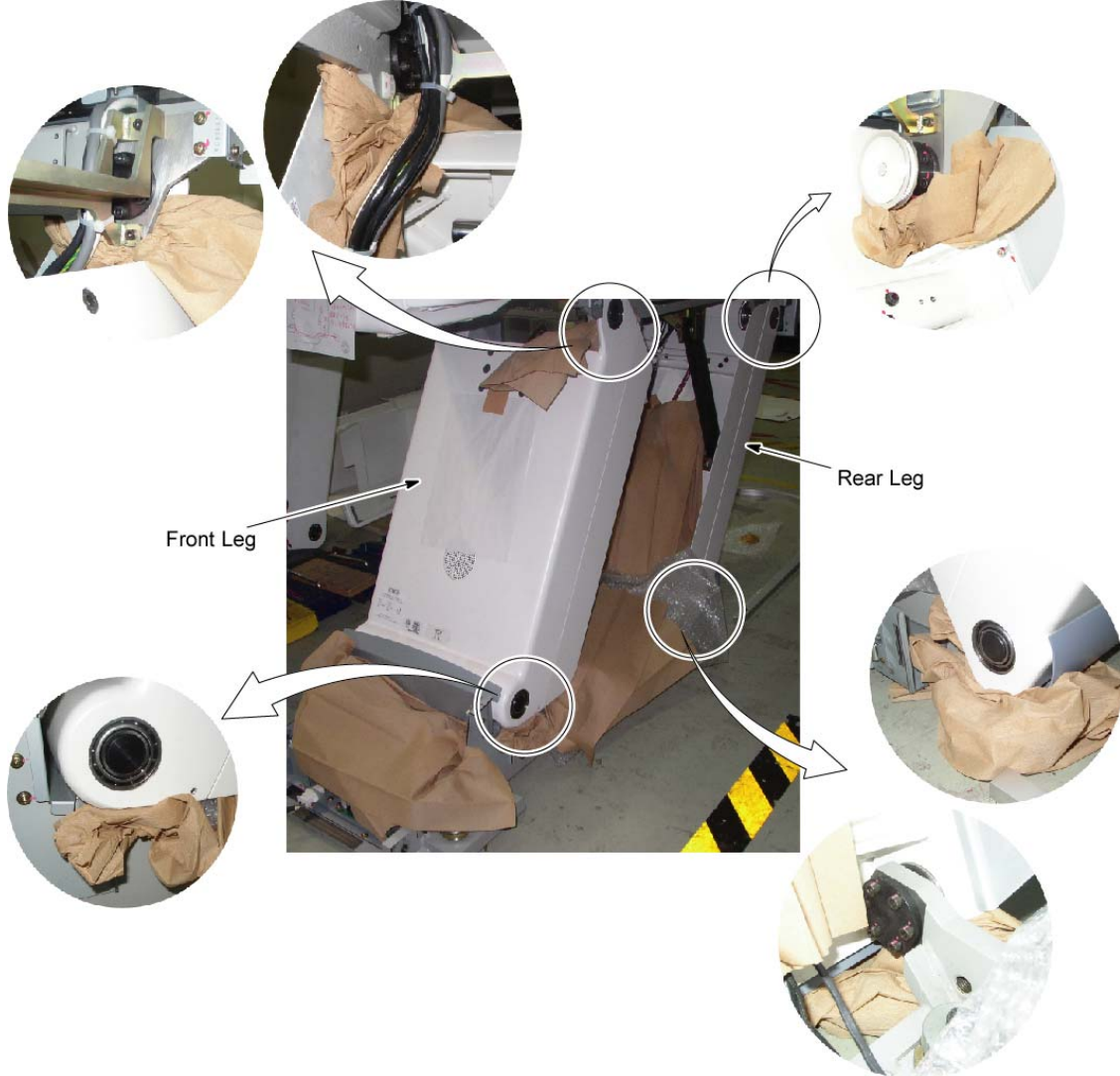


Figure - 6 Table Covering

GE Healthcare MI&CT Class C eService Notes are copyright (c) protected to General Electric Company Inc with all rights reserved. Class C eService Notes issued by GE Healthcare are proprietary to GE Healthcare Systems. Disclosure to persons not employed by GE Healthcare or not in possession of a Advanced Service Agreement with GE Healthcare is prohibited.

Duplication and distribution of any class of eService Notes, in any form by entities other than GE Healthcare is strictly prohibited.

- 10). Spray a leg bearing with CRC5-56 of equivalent spray lubricant (WD-40, etc.) to flush debris from each bearings, and wipe off the oil flowed from the bearing.

(If necessary, change the Kim-wipe or waste cloth and perform this step repeatedly.)

Spray Lubricant (CRC5-56, WD-40, etc.)



Figure - 7 Spraying

If apparent debris from the bearing appears, please write a PQR to report this problem.

Debris

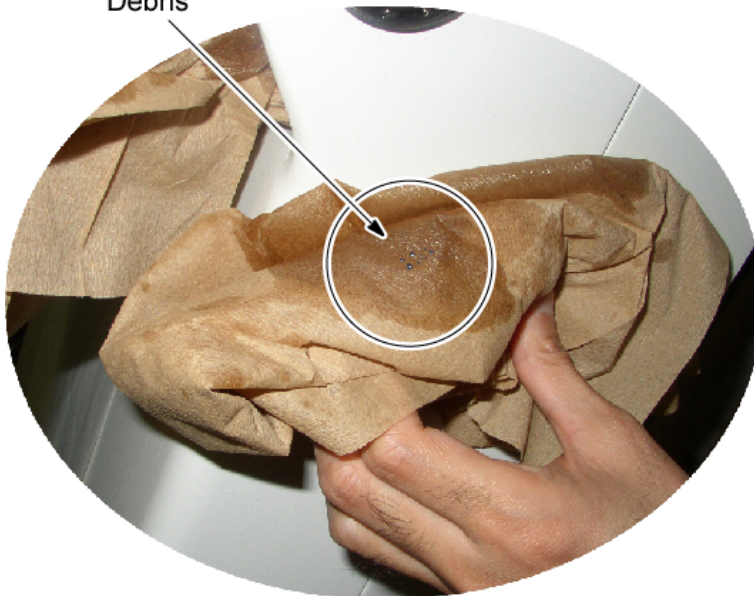


Figure - 8 Debris

Duplication and distribution of any class of eService Notes, in any form by entities other than GE Healthcare is strictly prohibited.

- 11). Visually check the retainer of the bearing for damage and deformation.

If the retainer is NOT OK, please report this problem by e-mail (fmi26812@ge.com).

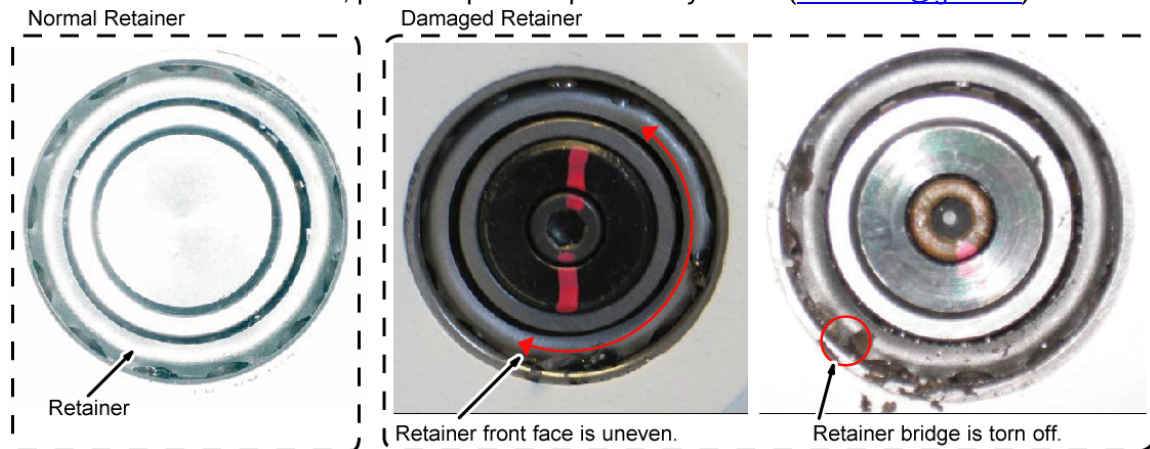


Figure - 9 Retainer

- 12). Flush other leg bearings as in steps 10 and 11.

- 13). Grease eight (8) leg bearings on both sides of the table.

Note: Firmly pack grease into the bearing so that it is applied to the inner race of the bearing, using a wooden stick or tie-wrap.

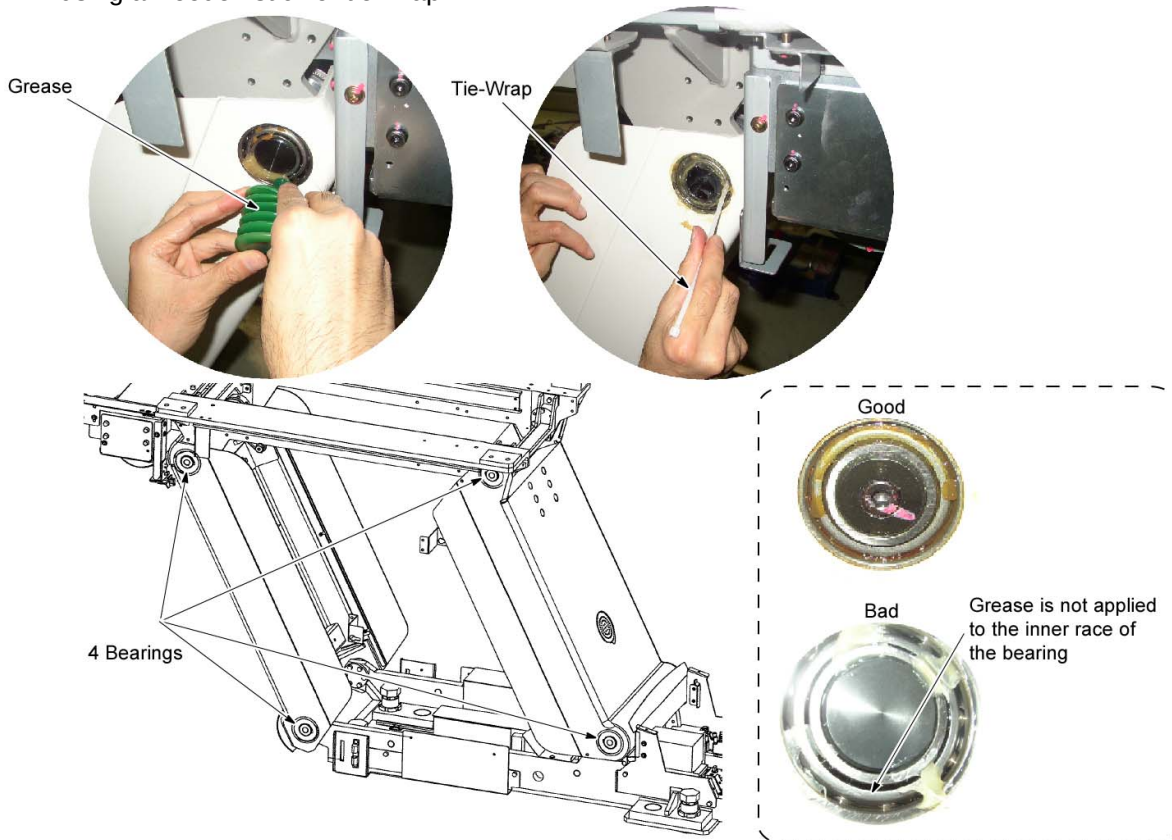


Figure - 10 Table Leg Bearing (Right Side)

Duplication and distribution of any class of eService Notes, in any form by entities other than GE Healthcare is strictly prohibited.

- 14). Move the table up and down several times to apply grease to bearings.
- 15). Wipe off any excess grease.
- 16). Power up the table using the table circuit breaker.
- 17). Verify that the table moves up and down normally.
- 18). Remove power from the table circuit breaker, and re-install cover plates and side plates as follows:
 - a.) Power up the table using the table circuit breaker, and lower the table to the down limit position, then remove power from the table circuit breaker.
 - b.) Attach the cover plate to the table top frame, and tighten two screws so that there is a clearance of 0.5 millimeter between the cover plate and rear leg.

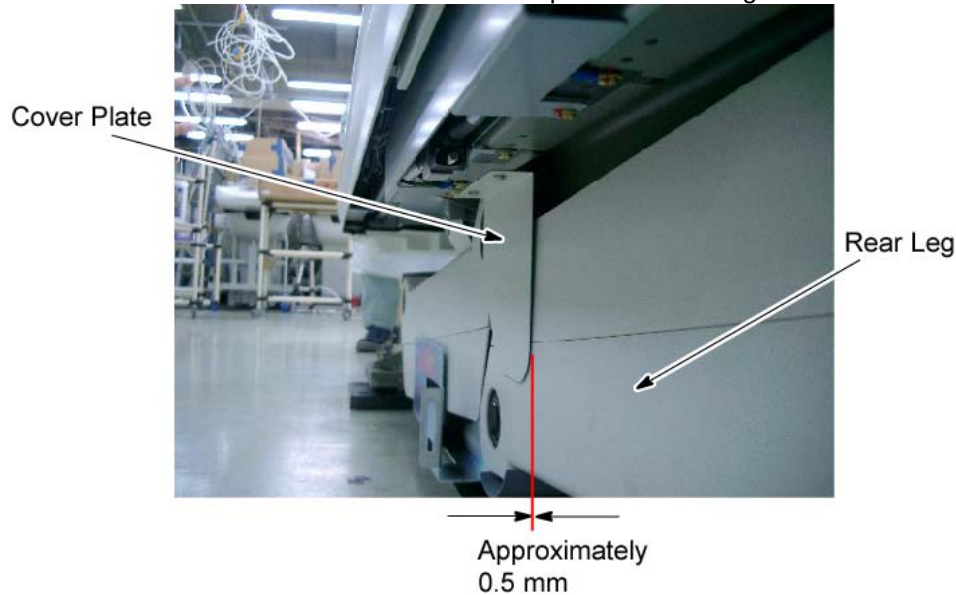


Figure - 11 Cover Plate Attachment

- c.) Power up the table using the table circuit breaker, and raise and lower the table repeatedly to verify that the cover plate does not interfere with the rear leg plate (Fig.A) and rear leg (Fig. B). Pay special attention to the points circled below.

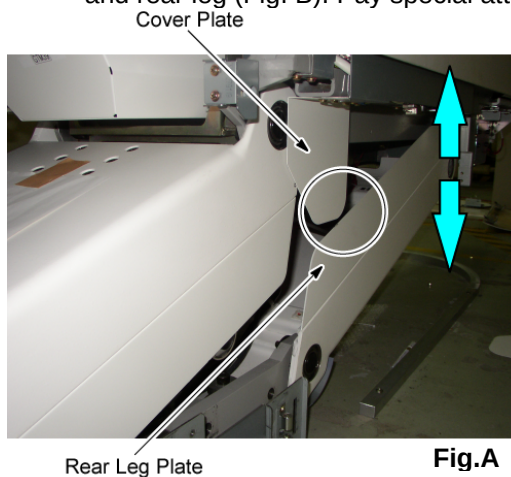


Fig.A

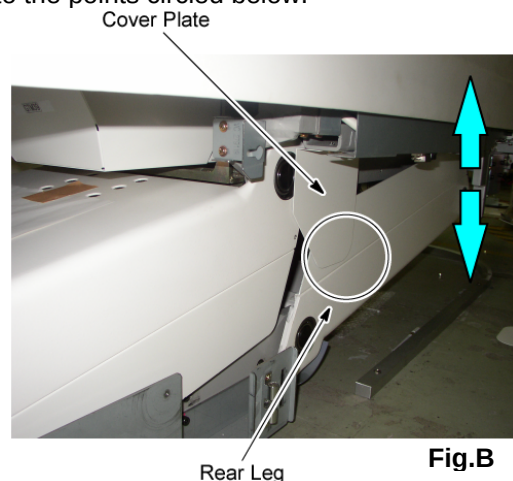


Fig.B

Figure - 12 Cover Plate and Rear Leg

Duplication and distribution of any class of eService Notes, in any form by entities other than GE Healthcare is strictly prohibited.

- d.) Raise the table to its highest position, and remove power from the table circuit breaker.
- e.) Attach the upper support bracket of the rear side plate to the table top frame, and hand-tighten the two screws.

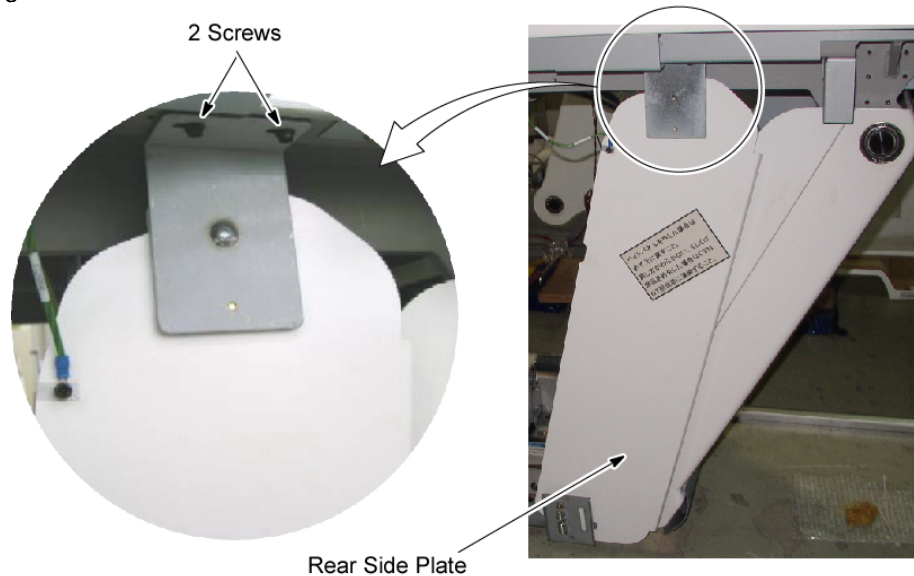


Figure - 13 Upper Support Bracket of the Rear Side Plate Attachment

- f.) Attach a spacer, a bush, and a shaft to the lower side of the rear side plate, and screw the side plate to the table base frame.

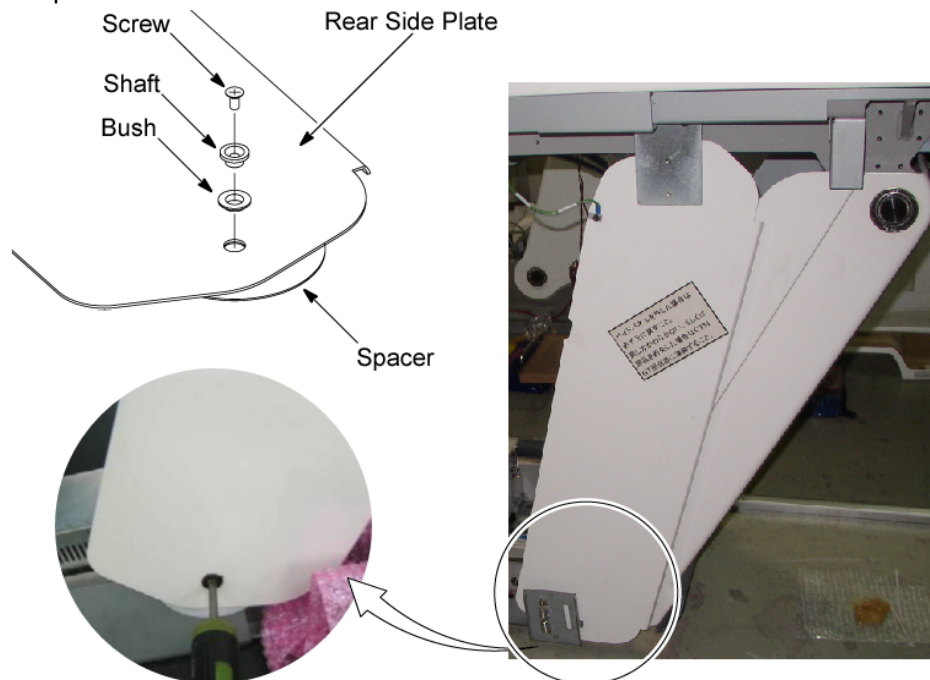


Figure - 14 Lower Side of the Rear Side Plate Attachment

Duplication and distribution of any class of eService Notes, in any form by entities other than GE Healthcare is strictly prohibited.

- g.) Push up the upper support bracket of the rear side plate against the table top frame, and pull the rear side plate away from the rear leg plate, then tighten the 2 screws of the upper support bracket of the rear side plate.

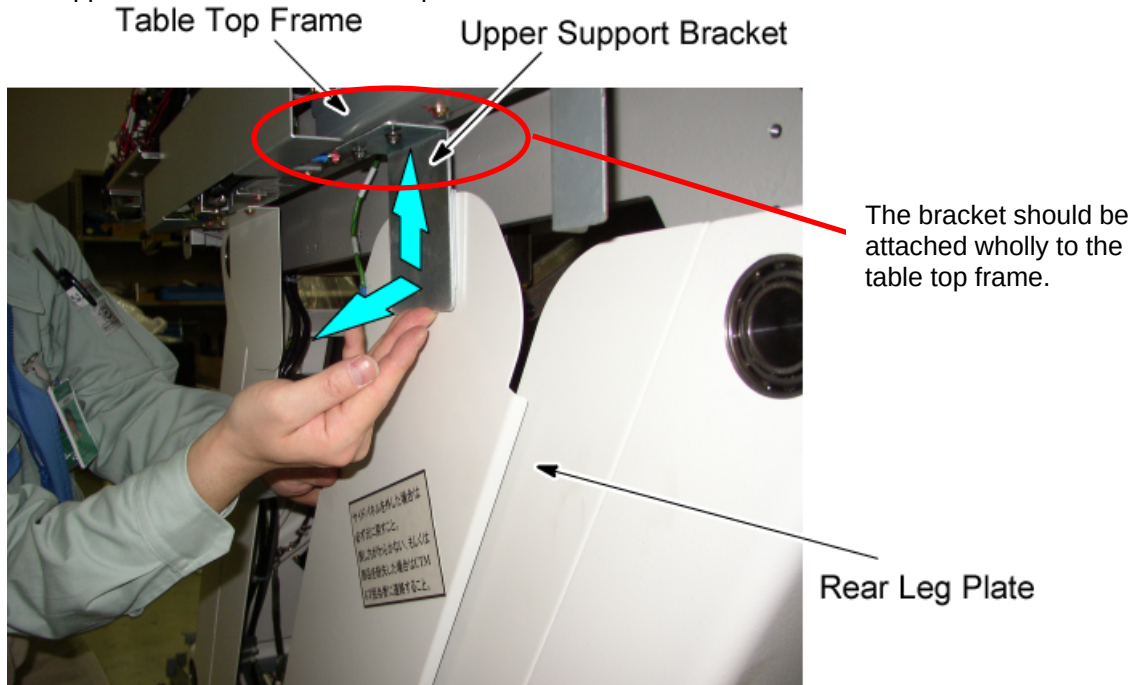


Figure - 15 Rear Side Plate Positioning

- h.) Connect the GND cable to the table frame.

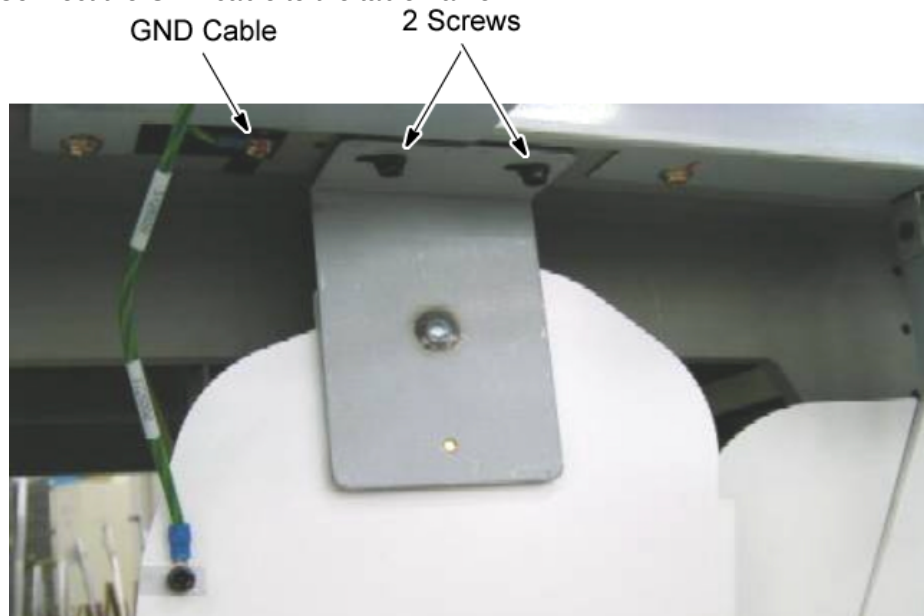


Figure - 16 Upper Support Bracket of the Rear Side Plate Fixing

Duplication and distribution of any class of eService Notes, in any form by entities other than GE Healthcare is strictly prohibited.

- i.) Power up the table using the table circuit breaker, and lower and raise the table repeatedly to verify that the rear side plate does not interfere with the leg plate.

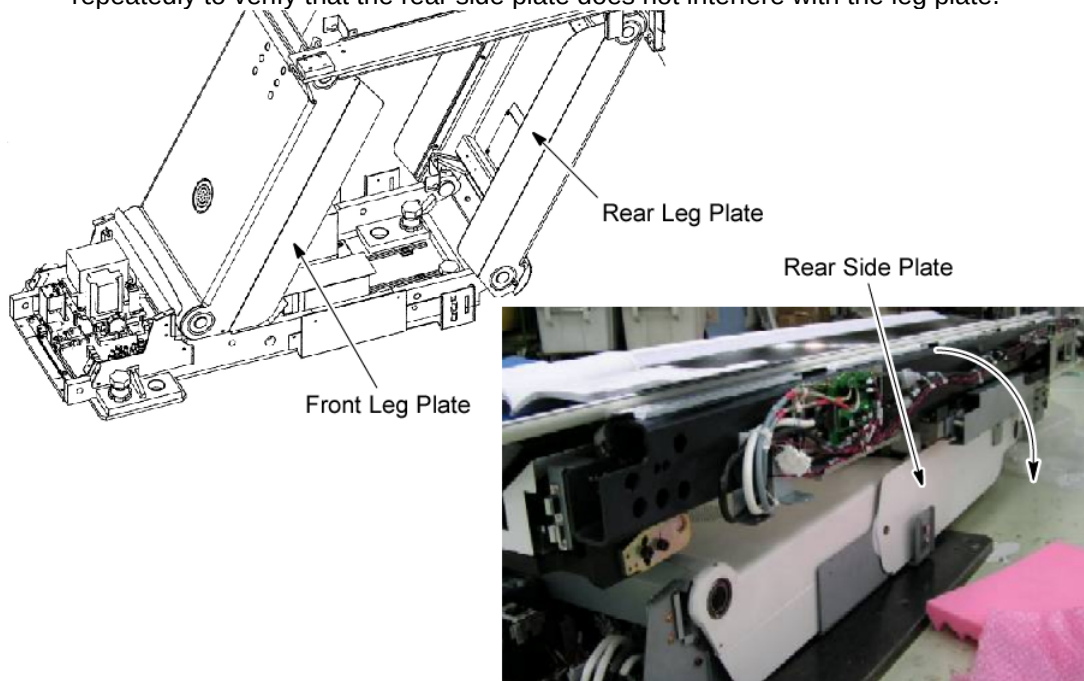


Figure - 17 Table Up/Down Check for Rear Side Plate

- j.) Raise the table to its highest position, and remove power from the table circuit breaker.
- k.) Attach the upper bracket of the front side plate to the table top frame, and hand-tighten the two screws.

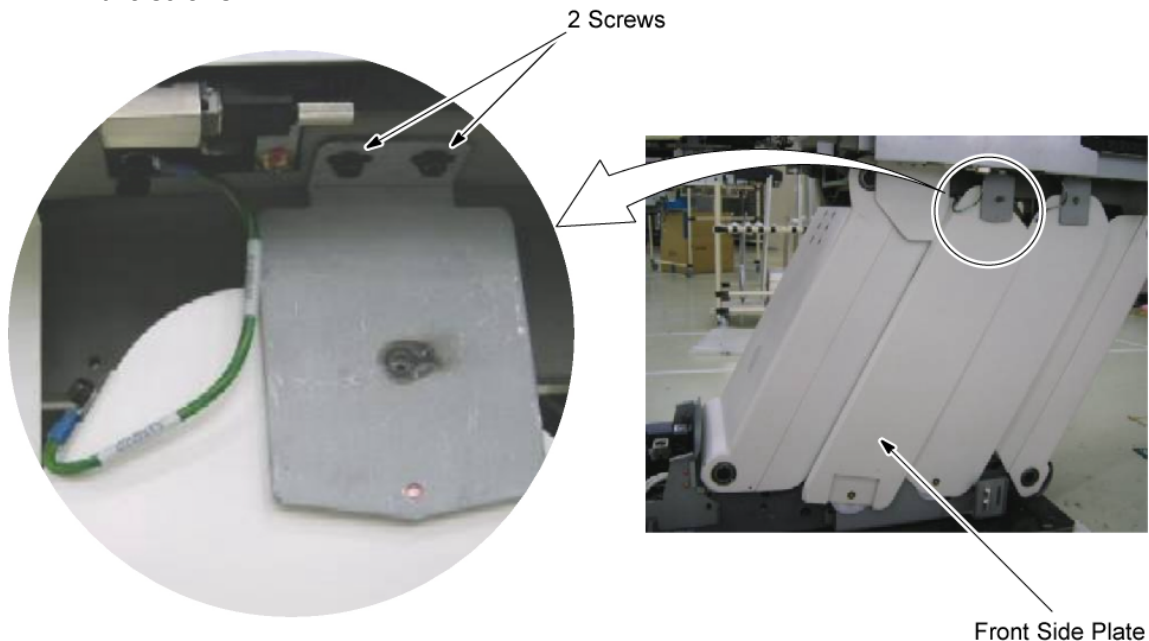


Figure - 18 Upper Support Bracket of the Front Side Plate Attachment

Duplication and distribution of any class of eService Notes, in any form by entities other than GE Healthcare is strictly prohibited.

- l.) Attach a spacer, a bush, and a shaft to the lower side of the front side plate, and screw the side plate to the table base frame

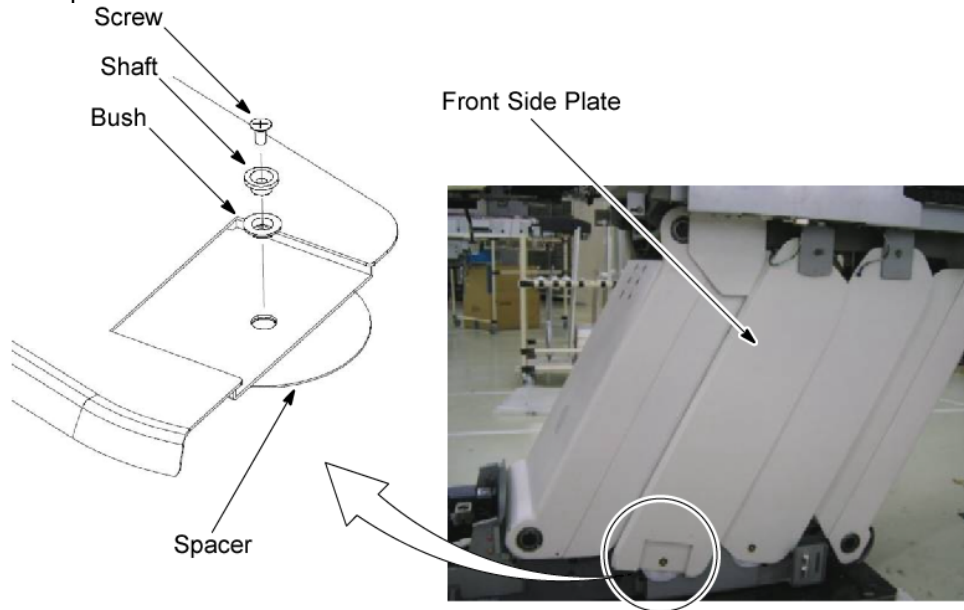


Figure - 19 Lower Side of the Front Side Plate Attachment

- m.) Push up the upper support bracket of the front side plate against the table top frame, and pull the upper support bracket away from the cover plate, then tighten the two (2) screws of the upper support bracket of the front side plate.
n.) Connect the GND cable to the table frame.

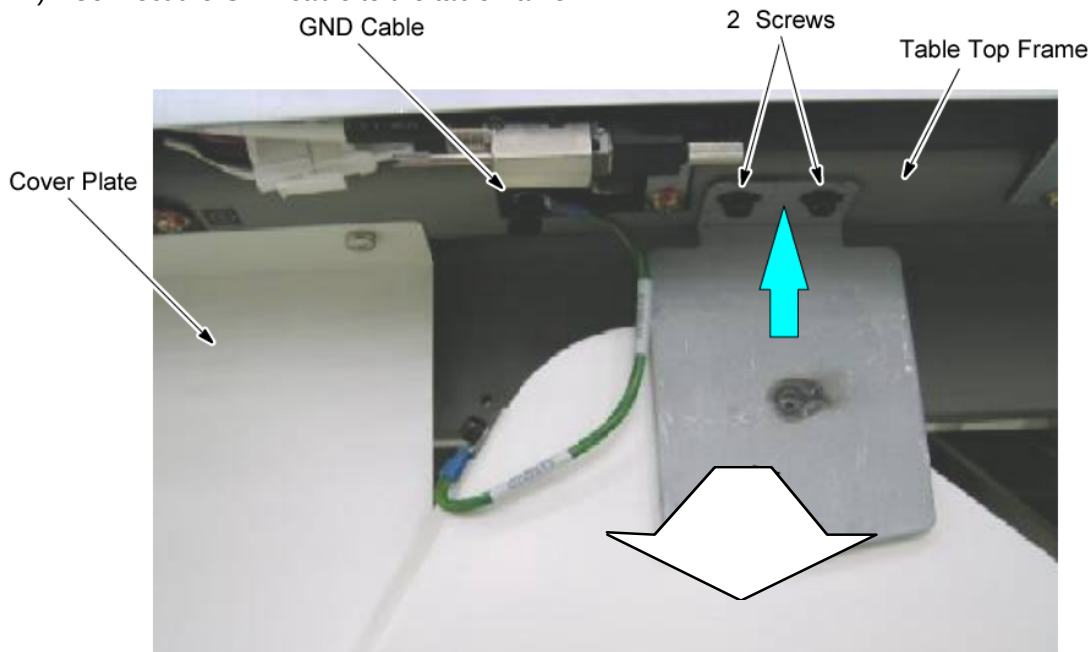


Figure - 20 Upper Support Bracket of the Front Side Plate Fixing

Duplication and distribution of any class of eService Notes, in any form by entities other than GE Healthcare is strictly prohibited.

- o.) Verify that the front side plate does not interfere with the cover plate.

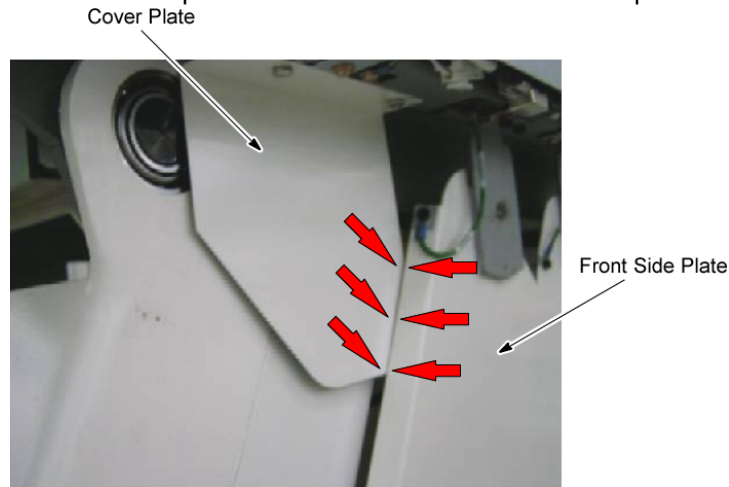


Figure - 21 Cover Plate and Front Side Plate

- p.) Power up the table using the table circuit breaker, and lower and raise the table repeatedly to verify that the front side plate does not interfere with the leg plate.

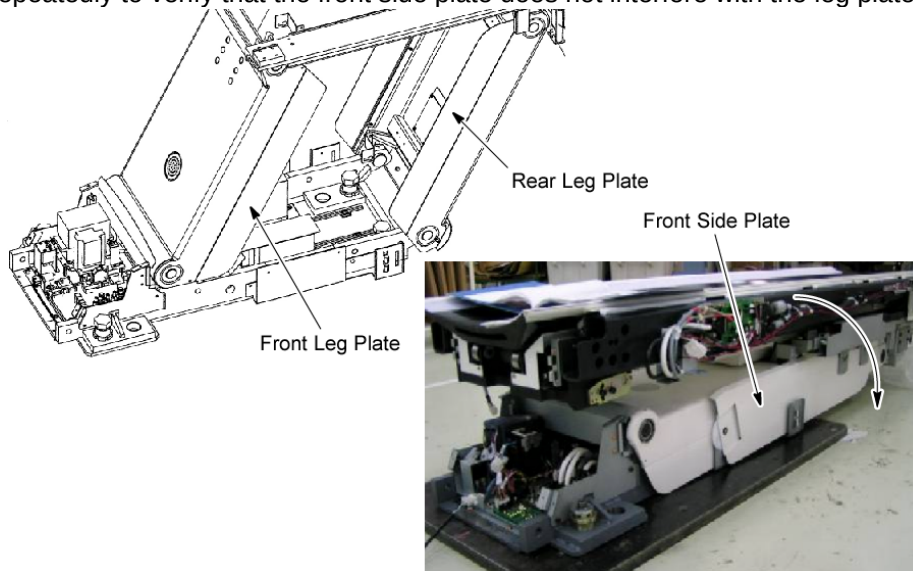


Figure - 22 Table Up/Down Check for Front Side Plate

- 19). Repeat step 18 for the other cover plate and side plates.
20). Remove power from the table circuit breaker, and re-install the front base cover, side base cover, IMS side covers, and IMS rear cover.

Takashi Uehara
Lead Service Integrator, GEYMS MI&CT Engineering

MI&CT Engineering Manager Signature: Mark Nyman