

GE Healthcare MI&CT Class C eService Notes are copyright (c) protected to General Electric Company Inc with all rights reserved. Class C eService Notes issued by GE Healthcare are proprietary to GE Healthcare Systems. Disclosure to persons not employed by GE Healthcare or not in possession of a Advanced Service Agreement with GE Healthcare is prohibited.
Duplication and distribution of any class of eService Notes, in any form by entities other than GE Healthcare is strictly prohibited.

**MI&CT ESERVICE NOTE
SN07062210a**

Class C information is Property of GE.

This short-term informational note is invalid 6 months after the issue date.

Date: June 22, 2007


Subject: Performix X-Ray Tube harness (2151264) may have loose wire crimps.

Effectivity: All Performix Tubes manufactured between April 2nd, 2007 and June 18, 2007 model numbers D3132T, D3139T, D3152T, D3176T, D3179T, D3185T, D3186T, D3189T.

Details:

There is a possibility that units listed in the effectivity above fail to produce x-rays due to a wiring failure in the harness 2151264. It is possible that the unit reaches the maximum allowed temperature and activates the protection circuit if the cooling system (oil pump and fans) is not energized properly. It is also possible, even though not reported to date, that the loose wires cause a short to the surrounding metal surfaces and cause the protection fuse to open in the CT gantry base.

Examples of the Harness defect

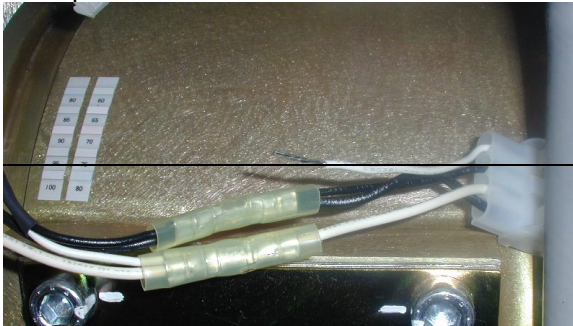

Resolution:

1. All units in the effectivity list above manufactured after June 18, 2007 have the harness 2151264 inspected for loose wires (pull test at supplier) and are ready for installation and use without inspection.
2. Inspect all units on hand, in effectivity list above, prior to installation.

GE Healthcare MI&CT Class C eService Notes are copyright (c) protected to General Electric Company Inc with all rights reserved. Class C eService Notes issued by GE Healthcare are proprietary to GE Healthcare Systems. Disclosure to persons not employed by GE Healthcare or not in possession of a Advanced Service Agreement with GE Healthcare is prohibited.

Duplication and distribution of any class of eService Notes, in any form by entities other than GE Healthcare is strictly prohibited.

Procedure

Harness 2151264 Manual Pull Test

Note: This procedure must be executed prior to X-Ray tub installation and electrical connection to the gantry. Do not use pliers to perform this procedure


1. Hold the AMP splicers in the middle and manually pull each of the 8 wires shown below (4 White & 4 Black).



2. If the test above fails, order a new tube unit and return the failed one as FOA.

Note: Do not attempt to repair the crimp in the field. The proper crimping tool for this splic connector is not available for field use.

Note: There is no need to pull-test the wires going to the 3 connectors shown below

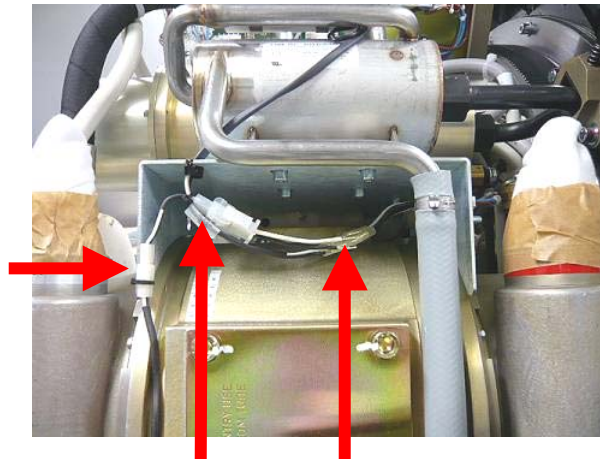

Affected Service Publications: Not identified by author. See Model numbers under “effectivity”.

PQR/SPR or other Tracking #:

Luiz Coutinho

GEHC/ISS Performix Product Champion

MI&CT Engineering Manager Reviewer: _____ Mark Nyman _____