

MI&CT eSERVICE NOTE

SN08100112a

Class C information is Property of GE.

Date: October 1, 2008

Subject: Keyboard and SCIM Deviation for GOC6 Operator Console

Effectivity: GOC6 Consoles (see below)

Console Part Numbers Affected

System Alias and Scanner Type	Console P/Ns affected
Discovery CT750HD	5212920-100
Lightspeed VCT 64	5212920-300

Details:

Due to Supply Chain issues, the new GOC6 Console will not have the black colored SCIM and Keyboards as per design. Instead the new console will be supplied with existing beige colored SCIM and Keyboards. Anticipated arrival of the black colored SCIM and Keyboards is 1st quarter 2009.

Consequently there is Service Documentation impact related to the use of the beige colored SCIM and Keyboards with the GOC6 Console. FRU part numbers and console keyboard installation are affected. *(Note: Beige colored SCIM installs the same as the black colored version on GOC6.)*

Plans are in place to upgrade the GOC6 Console installed base with black SCIM and Keyboards once available.

GE Healthcare MI&CT Class C eService Notes are copyright (c) protected to General Electric Company Inc with all rights reserved. Class C eService Notes issued by GE Healthcare are proprietary to GE Healthcare Systems. Disclosure to persons not employed by GE Healthcare or not in possession of a Advanced Service Agreement with GE Healthcare is prohibited.
Duplication and distribution of any class of eService Notes, in any form by entities other than GE Healthcare is strictly prohibited.

Resolution:

FRU Part Numbers:

Beige SCIM 2262158

Beige Keyboards


5271210-EN	English
5271210-DA	Danish
5271210-NL	Dutch
5271210-FR	French
5271210-DE	German
5271210-IT	Italian
5271210-NO	Norwegian
5271210-ES	Spanish
5271210-SV	Swedish
5271210-PT-PT	Portuguese, European
5271210-PT-BR	Portuguese, Brazil

Keyboard Installation and Replacement:

The beige colored keyboards utilize a PS2 interface, verses the newer USB interface. With design changes implemented in the GOC6 Console, support for the older PS2 keyboard connector is no longer present on the Upper Bulkhead Panel connector. Consequently, the beige colored keyboards must be connected directly to the Host PC's PS2 Keyboard Port inside the GOC6 console. When installing the beige colored keyboards, follow the connection process below. Do not follow the instructions in the GOC6 Service Manuals. These procedures are for the new USB connected black keyboards.

GE Healthcare MI&CT Class C eService Notes are copyright (c) protected to General Electric Company Inc with all rights reserved. Class C eService Notes issued by GE Healthcare are proprietary to GE Healthcare Systems. Disclosure to persons not employed by GE Healthcare or not in possession of a Advanced Service Agreement with GE Healthcare is prohibited.

Duplication and distribution of any class of eService Notes, in any form by entities other than GE Healthcare is strictly prohibited.

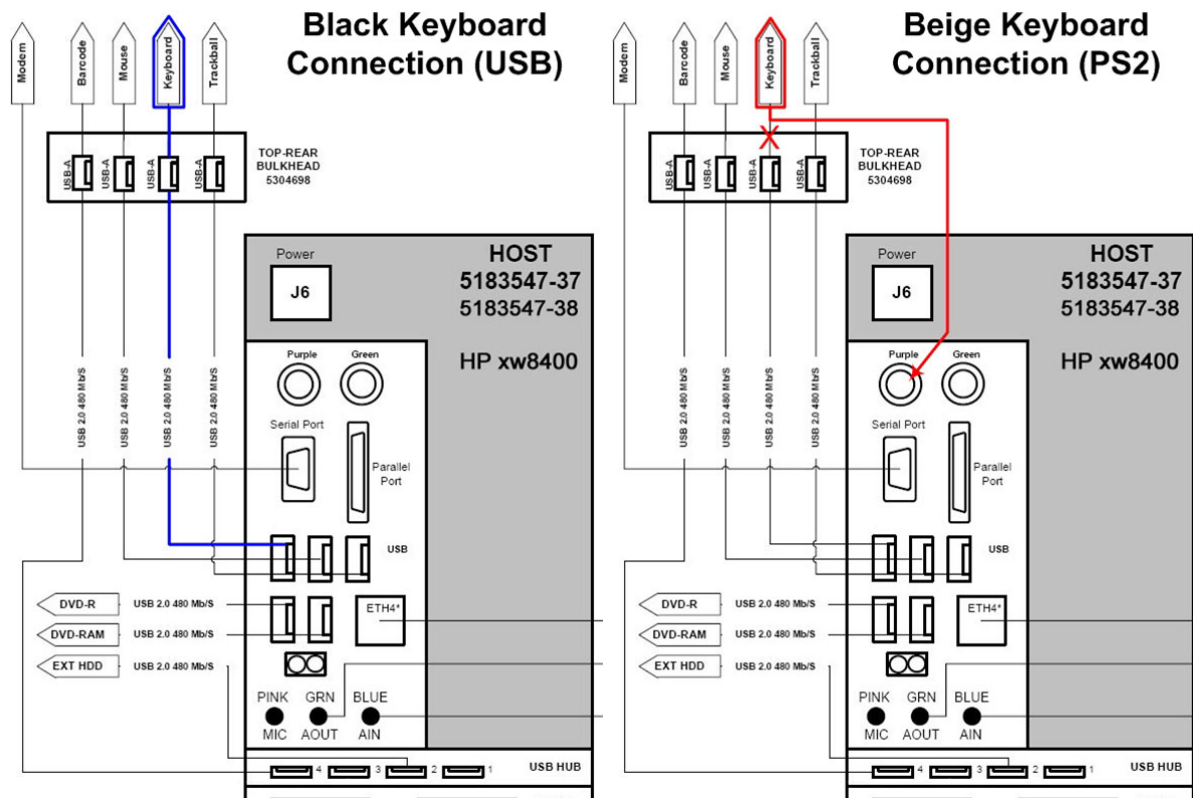

When using the beige colored PS2 style keyboards, connect the keyboard directly to the Host Computer's Purple PS2 Port on the rear of the HP xw8400 PC. See Illustration below.



Pass the keyboard cable along the Upper Rear Bulkhead panel at the top of the console cabinet and behind the console's rear cover assembly, along with the Monitor Video and DVD Peripheral Tower USB cables.

GE Healthcare MI&CT Class C eService Notes are copyright (c) protected to General Electric Company Inc with all rights reserved. Class C eService Notes issued by GE Healthcare are proprietary to GE Healthcare Systems. Disclosure to persons not employed by GE Healthcare or not in possession of a Advanced Service Agreement with GE Healthcare is prohibited.
Duplication and distribution of any class of eService Notes, in any form by entities other than GE Healthcare is strictly prohibited.

Affected Service Publication:

LightSpeed VCT64 w/GOC6 Supplemental	5318716-2EN, Revision 01
LightSpeed VCT64 w/GOC6 Supplemental (Advanced)	5321294-2EN, Revision 01
Discovery CT750HD - Service Information	5307450-2EN, Revision 01
Discovery CT750HD - Service Information (Advanced)	5307452-2EN, Revision 01

PQR/SPR #: None

Michael L. Bauer
Console Subsystem SMD
MI&CT Service Engineering