

[new...](#) | [new revision](#) | [page status](#) | [copy](#) | [move](#) | [delete](#)

welcome!
ct tech pubs web page
ct water cooler
ct tech pubs updates
surveys
9800/hla
highspeed family (us)
ct/i
lightspeed family
yms family
all ct products
index
old ct emails
communiques pending
approval
customize
members

search:

advanced search
news: daily | weekly

go offline | chat | notify |
print | tutorial | help

home> lightspeed family

Proper Signal Brush Placement

christopher hardiman (christopher c hardiman), 06/20/2002 - 06:35 am



PROPRIETARY WARNING

This information is Property of GE - For GE Service Personnel Only.
Disclosure to third parties prohibited.

Communique Number

MQUE-5B3P2R

[Revision #CH-A2FE of **Proper Signal Brush Placement by Christopher C Hardiman**]

Equipment

Schliefring Slipring on LightSpeed & HiSpeed Systems

Hints/Tips

The LightSpeed and HiSpeed Systems', using the Schliefring platter style slip ring, signal brush replacement kit (P/N: 2245480) contains two different types of signal brushes. Six of these (P/N: 2253678) have red-marked contact surfaces (connection end of the assembly) and contain a higher silver content for cleaning the signal rings. It is VITAL that these brushes be placed in the "Brush Type 1" Red Zone on the brush block during replacement (see 2240267ADW_r1.pdf, below). The location and higher silver content optimize the self-cleaning function and prolong the life of all of the signal brushes. Refer to the attached diagrams showing the differences in the brushes.

You can download any of the files below by dragging them to your Windows desktop.
To open a file for reading, double-click it.







Author

Christopher C Hardiman

Revision


06/20/2002 06:35 AM

[new...](#) | [new revision](#) | [page status](#) | [copy](#) | [move](#) | [delete](#)

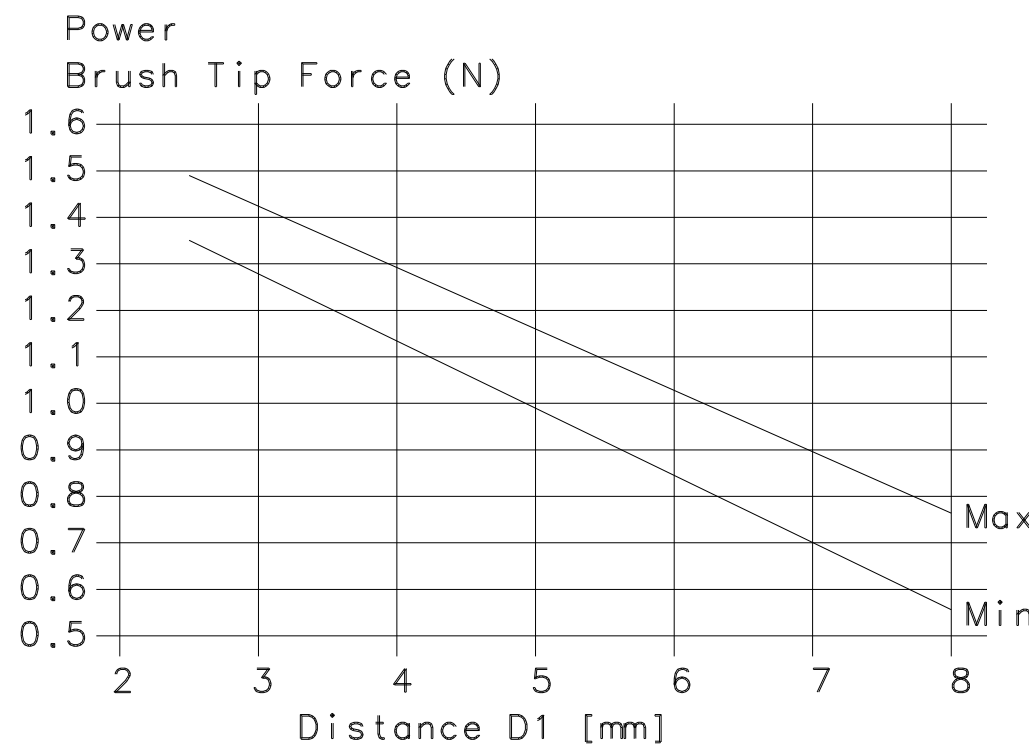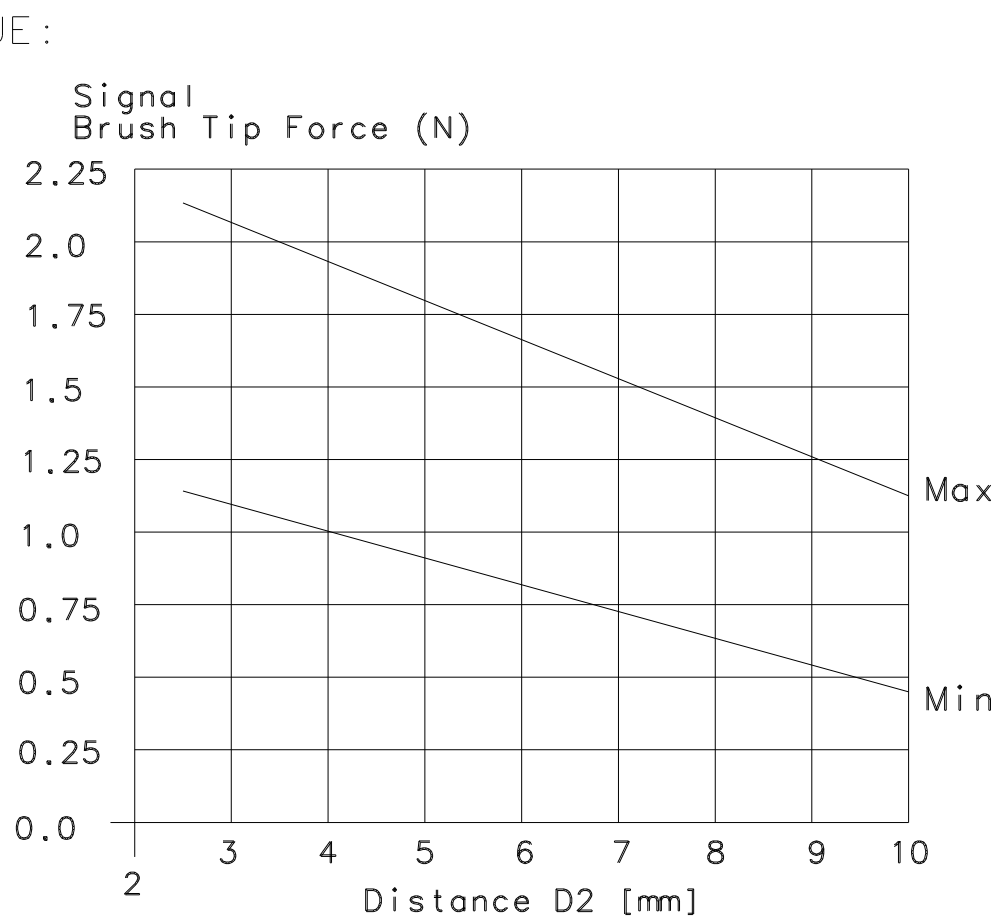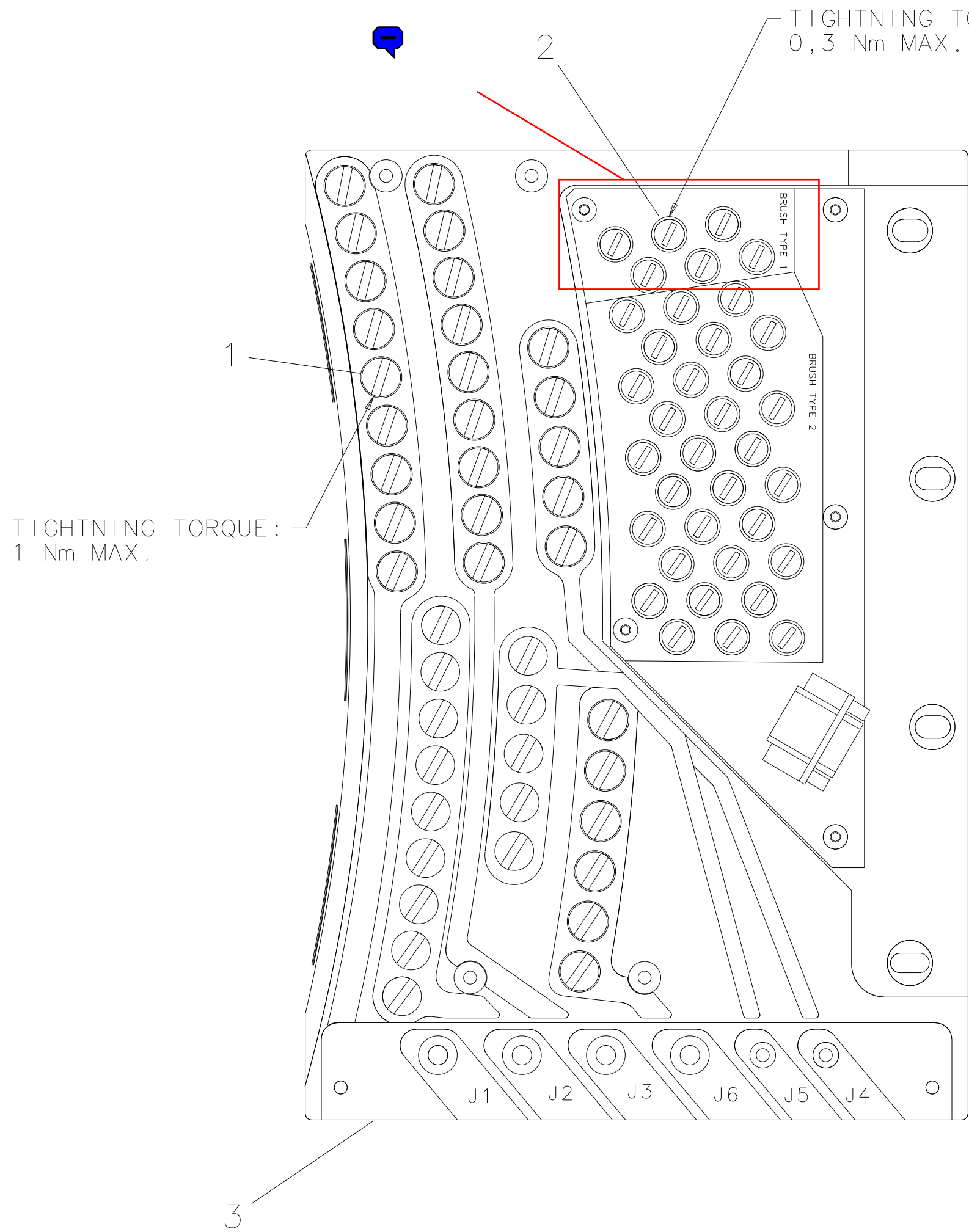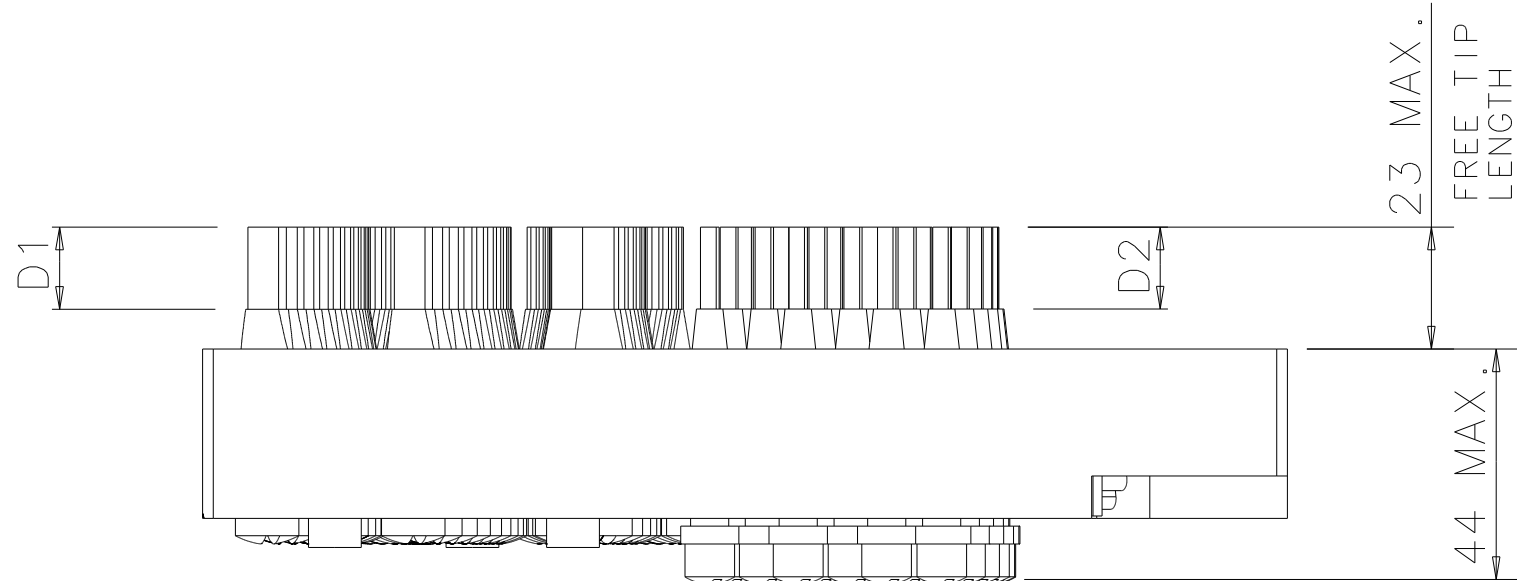


TIGHTNING TORQUE
M4 = 1,8 Nm MAX.
M6 = 6,0 Nm MAX.

REVISIONS						
REV	QTY	MADE BY	CH	DATE	PCN/FZ2 NO.	DESCRIPTION
B		RANDY WENDT		99SEP24	NCN	
REMOVED PIECE PART DIMENSIONS						
C		RANDY WENDT		00JAN24	NCN	
REVISED PER MOLDED DESIGN, ADD TIP DIMENSIONS AND TORQUES						
D		RANDY WENDT		00FEB11	NCN	
REVERSED ITEMS 1 & 2, SHOW TERM COVER, ADD (N) TO GRAPHS.						
1		JAIME LEHRER		02MAY13	220500	
UPDATE POWER BRUSH TIP FORCE GRAPH						

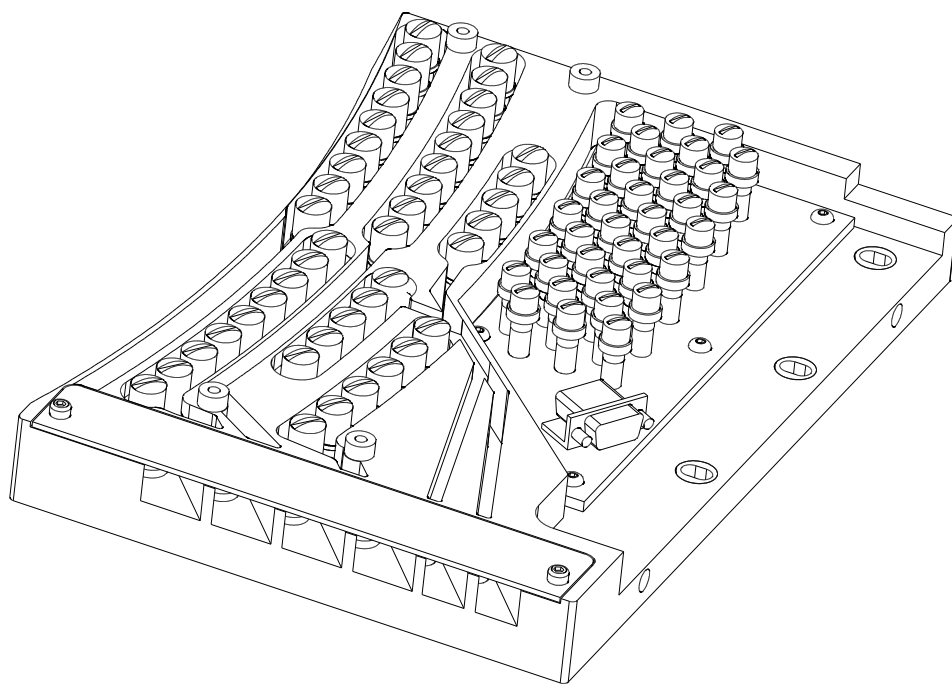

PHYSICAL ID: 2240267

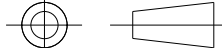MATERIAL:		FINISH:		SURFACE ROUGHNESS ✓		
ELECTRONIC SIGNOFF PRODUCT RELEASED		METRIC	APPLIED PRACTICES	THIRD ANGLE PROJECTION 		
				TITLE BRUSHBLOCK ASSY		
MADE BY RANDY WENDT		99APR13	GE MEDICAL SYSTEMS MILWAUKEE, WI		2240267ADW	
ISSUED E J CASTILLO		00FEB22			SHEET 1 OF 1	


REVISIONS

REV	QTY	MADE BY	CH	DATE	PCN/FZ2 NO.	DESCRIPTION
B		RANDY WENDT		00JAN06	NCN	

ADD WEAR MARK, ADD MATERIAL

C		RANDY WENDT		00FEB11	NCN	
---	--	-------------	--	---------	-----	--

4X3 WAS 3SQ., ADD YYWW MARKING

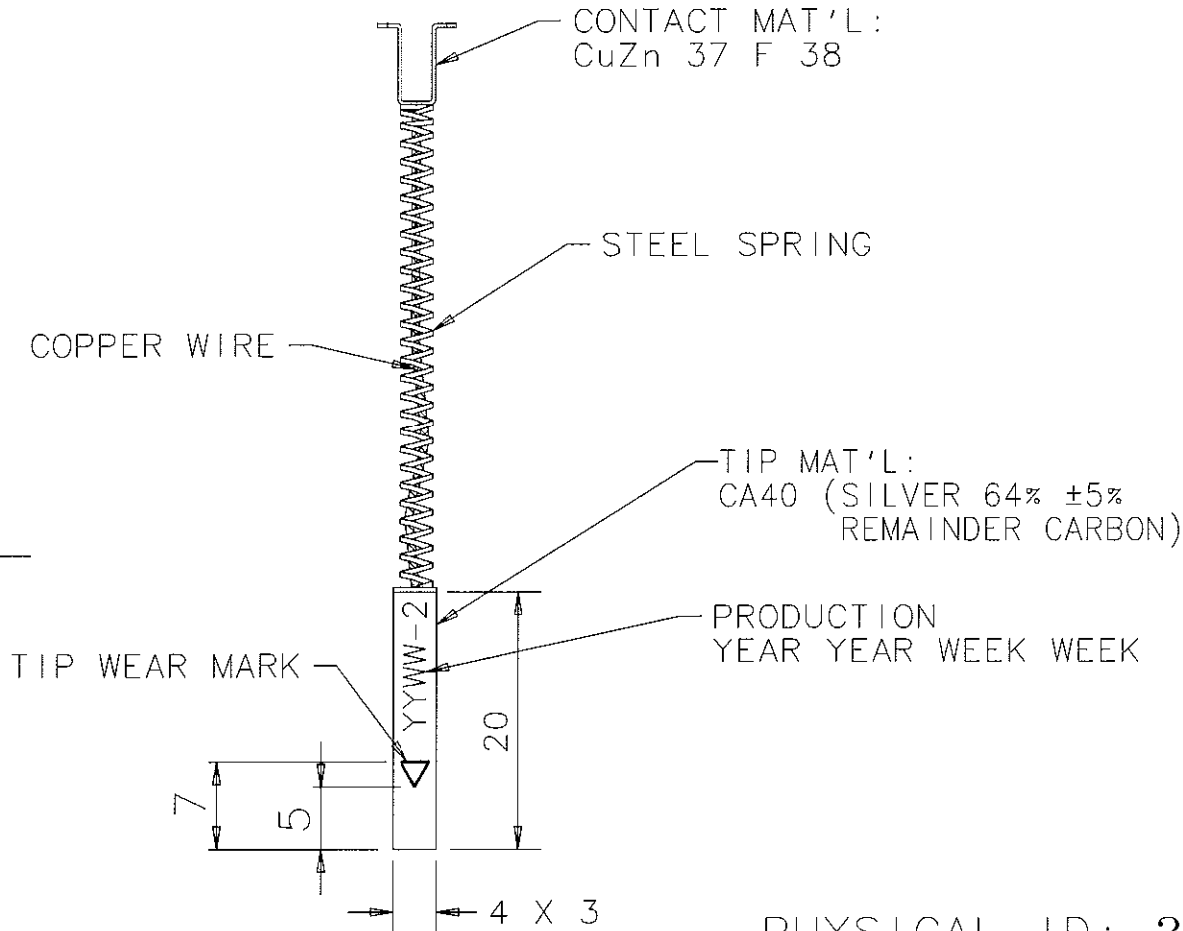

PHYSICAL ID: 2253677

MATERIAL:	FINISH:	SURFACE ROUGHNESS ✓
-----------	---------	------------------------

		METRIC	APPLIED PRACTICES 46-316917	THIRD ANGLE PROJECTION
CHECKER				
SOURCING				
MFG	JVL 00FEB15			
ENG	EJC 00FEB15			
DES	RNW 00FEB15			
APPROVALS	YY/MM/DD	TITLE SIGNAL BRUSHTIP		REV 0
MADE BY RANDY WENDT	99AUG19	FIRST MADE FOR: HELI-GAN-SLIPR		MOD CODE 2253677DDW
ISSUED/CH 	100221	GE MEDICAL SYSTEMS MILWAUKEE, WI		SHEET 1 OF 1

FIRST DESIGNED IN CONJUNCTION WITH:

REV C 98/JAN/03

SDRC
A4


REVISIONS

REV	QTY	MADE BY	CH	DATE	PCN/FZ2 NO.	DESCRIPTION
B		RANDY WENDT		00JAN06	NCN	

ADD TIP WEAR MARK , ADD RED MARKING, ADD MATERIALS



C		RANDY WENDT		00FEB11	NCN
---	--	-------------	--	---------	-----

4X3 WAS 3 SQ., ADD YYWW MARKING

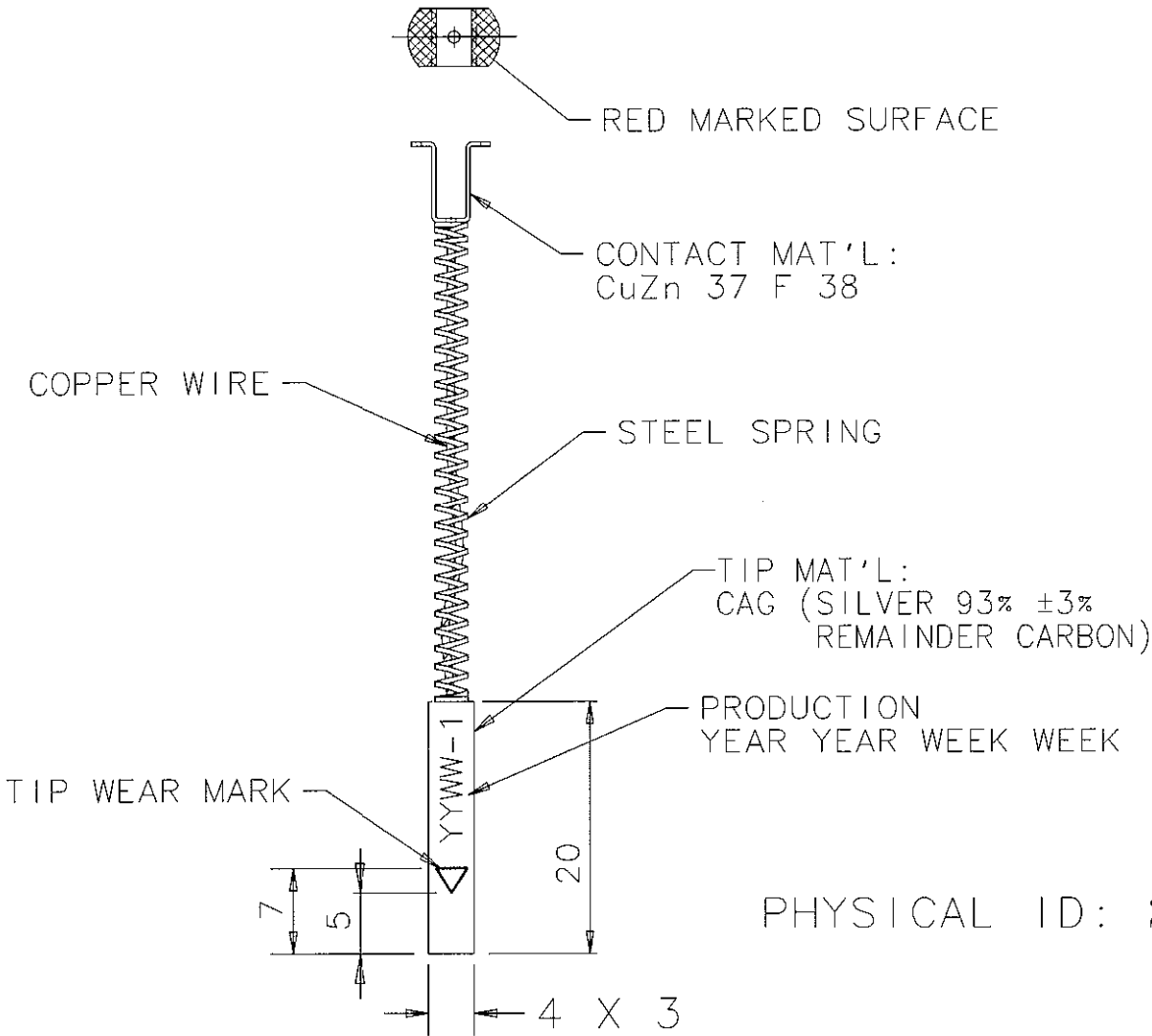

PHYSICAL ID: 2253678

MATERIAL:	FINISH:	SURFACE ROUGHNESS ✓
-----------	---------	------------------------

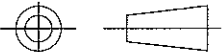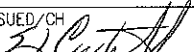		METRIC	APPLIED PRACTICES	THIRD ANGLE PROJECTION
CHECKER			46-316917	
SOURCING				
MFG	KVL 00FEB15			
ENG	EJC 00FEB15			
DES	RNW 00FEB15			
APPROVALS	YY/MM/DD	TITLE		REV
		SIGNAL BRUSHTIP		0
		FIRST MADE FOR: HELI-GAN-SLIPR		MOD CODE
MADE BY	RANDY WENDT 99AUG19	GE MEDICAL SYSTEMS		2253678DDW
ISSUED/CH	 000222	MILWAUKEE, WI		SHEET 1 OF 1

FIRST DESIGNED IN CONJUNCTION WITH:

REV C 96/JAN/03

SDRC
A4