



PATIENT HANDLING

TABLE OF CONTENTS

<u>SECTION</u>	<u>PAGE</u>
3 – PATIENT HANDLING	
Patient Handling Table Connections.....;	3-1
SRI4 Board Interconnects.....;	3-2
Longitudinal Drive Motor.....	3-3
Patient Blower Box.....	3-4

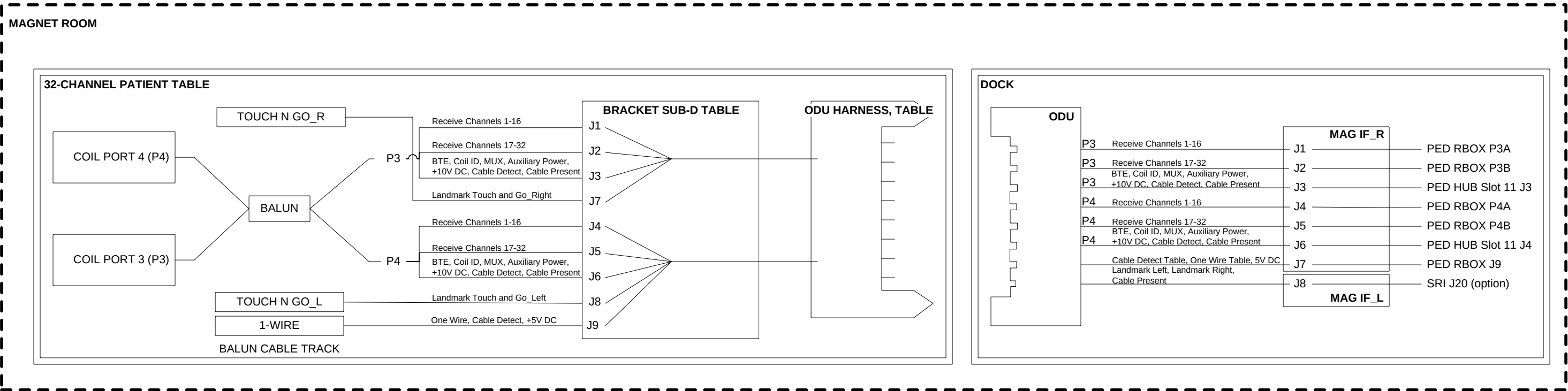


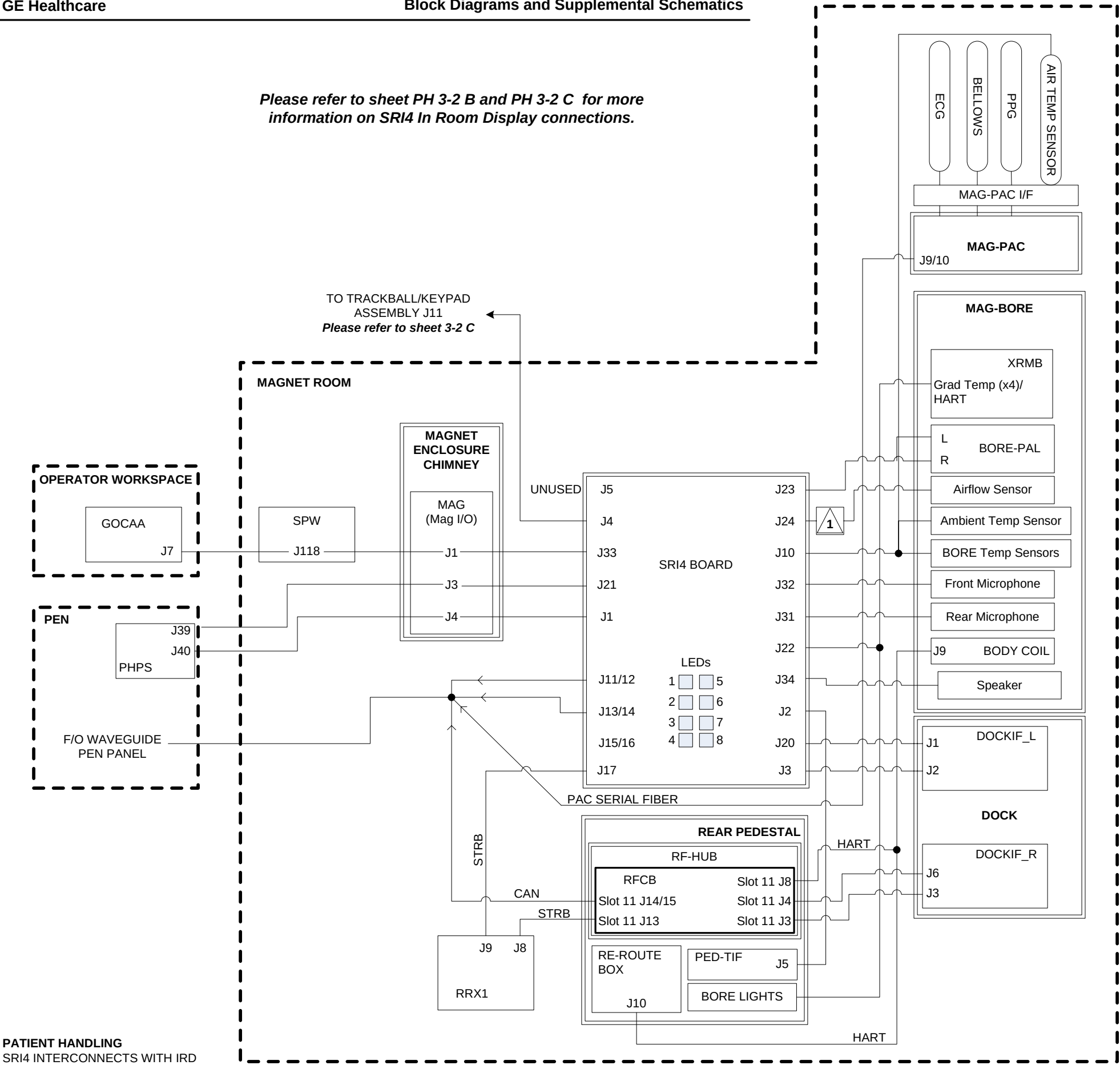

SRI CONNECTORS	CONNECTOR LABELS	SIGNALS
J1	SRI4 POWER (FROM PHPS)	VENT_CTRL_OUT 1W_INTERFACE ESTOP_BUS
J2	TABLE INTERFACE BOX	TABLE_FDBK 1W_DATA_INTF_BOX
J3	DOCK	DOCK_DN_CTRL DOCK_STAT
J4	CONTROL PANEL	PUSH_BUTTON_8V PAN_IND_8V DIM_LT_LED_8V ACTIVE_LED_8V BORE_LT_LED_8V BORE_VENT_LED_8V ESTOP_BUS ESTOP_PANEL_8V 1W_CTRL_PAN
J5	SEVEN SEGMENT DISPLAY	LONG_SSD_14V STR_SSD_14V LED_SSD_14V ESTOP_BUS ESTOP_DISPLAY
J10	BORE AND AMBIENT AIR TEMP	BORE_TEMP_SENSOR SCAN_ROOM_TEMP_SENS BORE_TEMP_CBL_PRES_N
J11	SRI CAN RX	FO_SRI_CAN
J12	SRI CAN TX	
J13	MCPDB CAN RX	FO_MCPDB_CAN
J14	MCPDB CAN TX	
J15	DEBUG TX	FO_DEBUG
J16	DEBUG RX	
J17	DABOUT_RX	FO_DABOUT
J20	LANDMARK (T AND GO)	LNDMRK_RIGHT LNDMRK_LEFT
J21	MCPDB (IN PHPS)	MCPDB_OUT MCPDB_STAT
J22	GRADIENT TEMP/BORE LIGHTS	XRBM_TEMP_SENSORS BORE_LIGHT_OUT XRBM_GRAD_TEMP_CBL_PRES_N 1W_DATA_XRBM_GRAD_TEMP BORE_LT_CBL_PRES_N
J23	PATIENT ALIGN LIGHTS (PAL)	PAL_CTRL_OUT PAL_CBL_PRES_N
J24	BODY COIL AIR FLOW SENSOR	AF_SENSORS AF_SENSORS_IF AF_CBL_PRES_N
J31	REAR MICROPHONE	RMIC
J32	FRONT MICROPHONE	FMIC
J33	CONSOLE INTERCOM	MIC_SIG SPK_SIG
J34	SPEAKER	SPK_SIG

SRI LED-NUMBERS	DESCRIPTION	FUNCTION
1	POWER OK	ON WHEN POWER IS ON
2	NIOS HEARTBEAT	FLASHES AT VARIABLE FREQUENCY
3	1W DIAG STATUS	ON WHEN SECONDARY DIAGNOSTIC NETWORK IS OK
4	ESTOP LOOP OPEN	ON WHEN ESTOP IS OPEN
5	MCPDB CAN STAT	FLASHES DURING NORMAL OPERATION
6	MCPDB CAN ERR	LATCHES AFTER CAN ERROR
7	SRI CAN STAT	FLASHES DURING NORMAL OPERATION
8	SRI CAN ERR	LATCHES AFTER CAN ERROR


1 Airflow Sensor Transducer Box

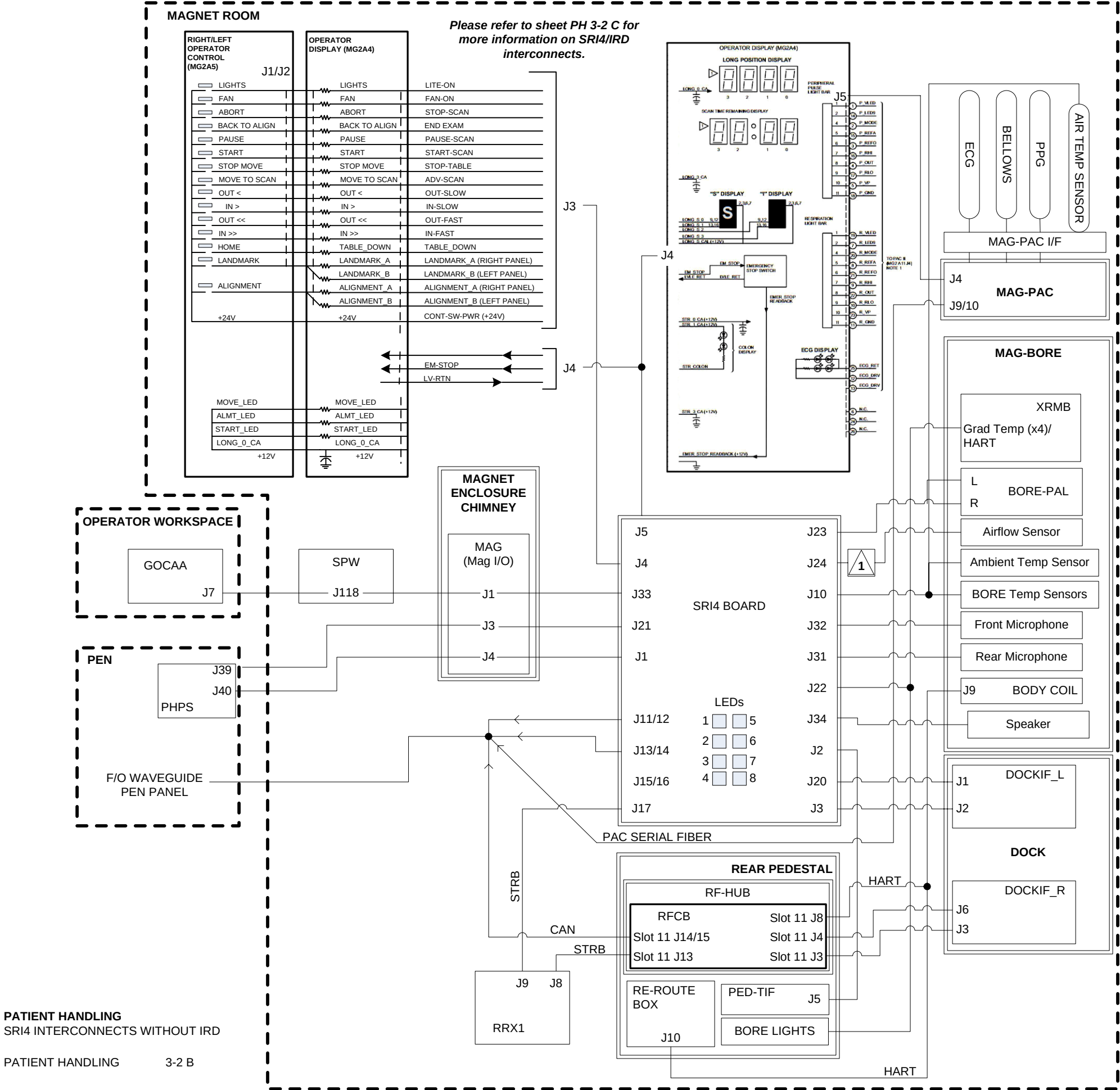


SRI CONNECTORS	CONNECTOR LABELS	SIGNALS
J1	SRI4 POWER (FROM PHPS)	VENT_CTRL_OUT 1W_INTERFACE ESTOP_BUS
J2	TABLE INTERFACE BOX	TABLE_FDBK 1W_DATA_INTF_BOX
J3	DOCK	DOCK_DN_CTRL DOCK_STAT
J4	CONTROL PANEL	PUSH_BUTTON_8V PAN_IND_8V DIM_LT_LED_8V ACTIVE_LED_8V BORE_LT_LED_8V BORE_VENT_LED_8V ESTOP_BUS ESTOP_PANEL_8V 1W_CTRL_PAN
J5	SEVEN SEGMENT DISPLAY	LONG_SSD_14V STR_SSD_14V LED_SSD_14V ESTOP_BUS ESTOP_DISPLAY
J10	BORE AND AMBIENT AIR TEMP	BORE_TEMP_SENSOR SCAN_ROOM_TEMP_SENS BORE_TEMP_CBL_PRES_N
J11	SRI CAN RX	FO_SRI_CAN
J12	SRI CAN TX	
J13	MCPDB CAN RX	FO_MCPDB_CAN
J14	MCPDB CAN TX	
J15	DEBUG TX	FO_DEBUG
J16	DEBUG RX	
J17	DABOUT RX	FO_DABOUT
J20	LANDMARK (T AND GO)	LNDMRK_RIGHT LNDMRK_LEFT
J21	MCPDB (IN PHPS)	MCPDB_OUT MCPDB_STAT
J22	GRADIENT TEMP/BORE LIGHTS	XRBM_TEMP_SENSORS BORE_LIGHT_OUT XRBM_GRAD_TEMP_CBL_PRES_N 1W_DATA_XRMB_GRAD_TEMP BORE_LT_CBL_PRES_N
J23	PATIENT ALIGN LIGHTS (PAL)	PAL_CTRL_OUT PAL_CBL_PRES_N
J24	BODY COIL AIR FLOW SENSOR	AF_SENSORS AF_SENSORS_IF AF_CBL_PRES_N
J31	REAR MICROPHONE	RMIC
J32	FRONT MICROPHONE	FMIC
J33	CONSOLE INTERCOM	MIC_SIG SPK_SIG
J34	SPEAKER	SPK_SIG

SRI LED-NUMBERS	DESCRIPTION	FUNCTION
1	POWER OK	ON WHEN POWER IS ON
2	NIOS HEARTBEAT	FLASHES AT VARIABLE FREQUENCY
3	1W DIAG STATUS	ON WHEN SECONDARY DIAGNOSTIC NETWORK IS OK
4	ESTOP LOOP OPEN	ON WHEN ESTOP IS OPEN
5	MCPDB CAN STAT	FLASHES DURING NORMAL OPERATION
6	MCPDB CAN ERR	LATCHES AFTER CAN ERROR
7	SRI CAN STAT	FLASHES DURING NORMAL OPERATION
8	SRI CAN ERR	LATCHES AFTER CAN ERROR

1 Airflow Sensor Transducer Box

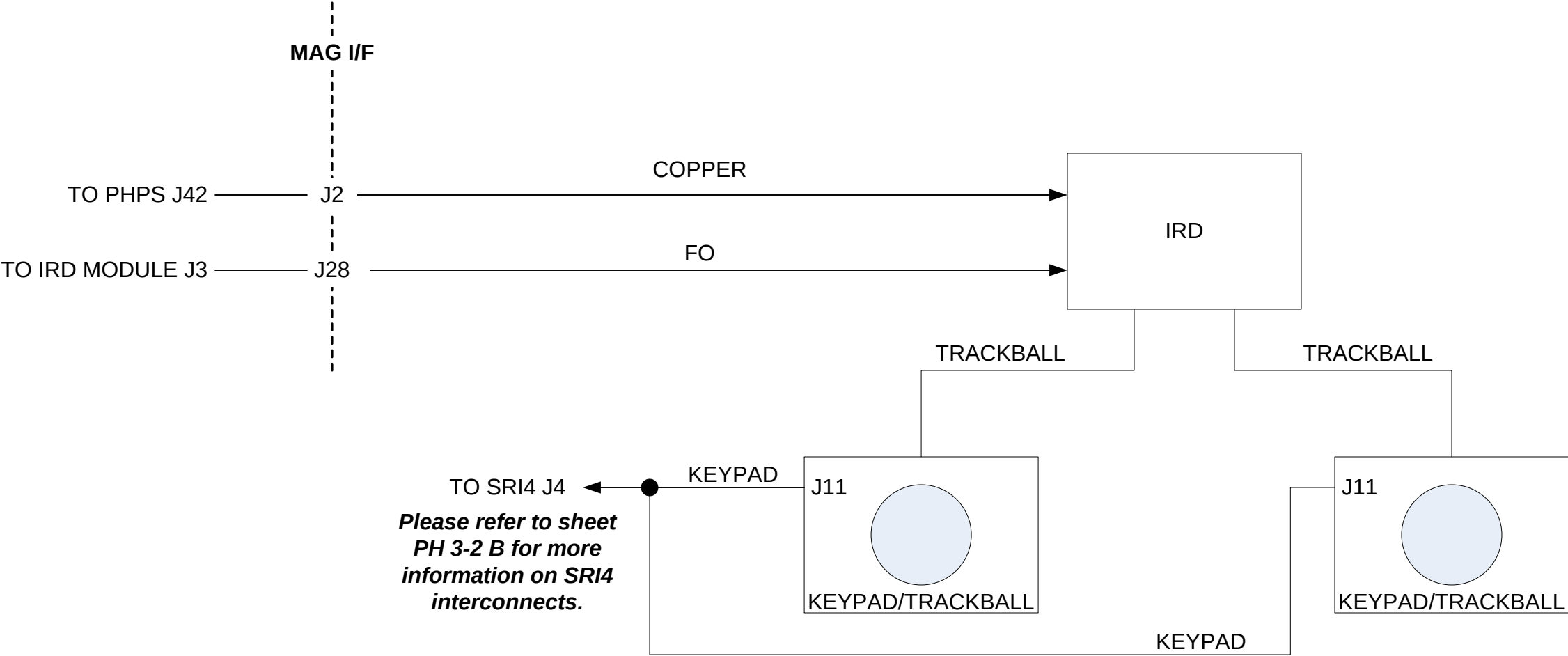

PATIENT HANDLING
In Room Display Connect to SRI4

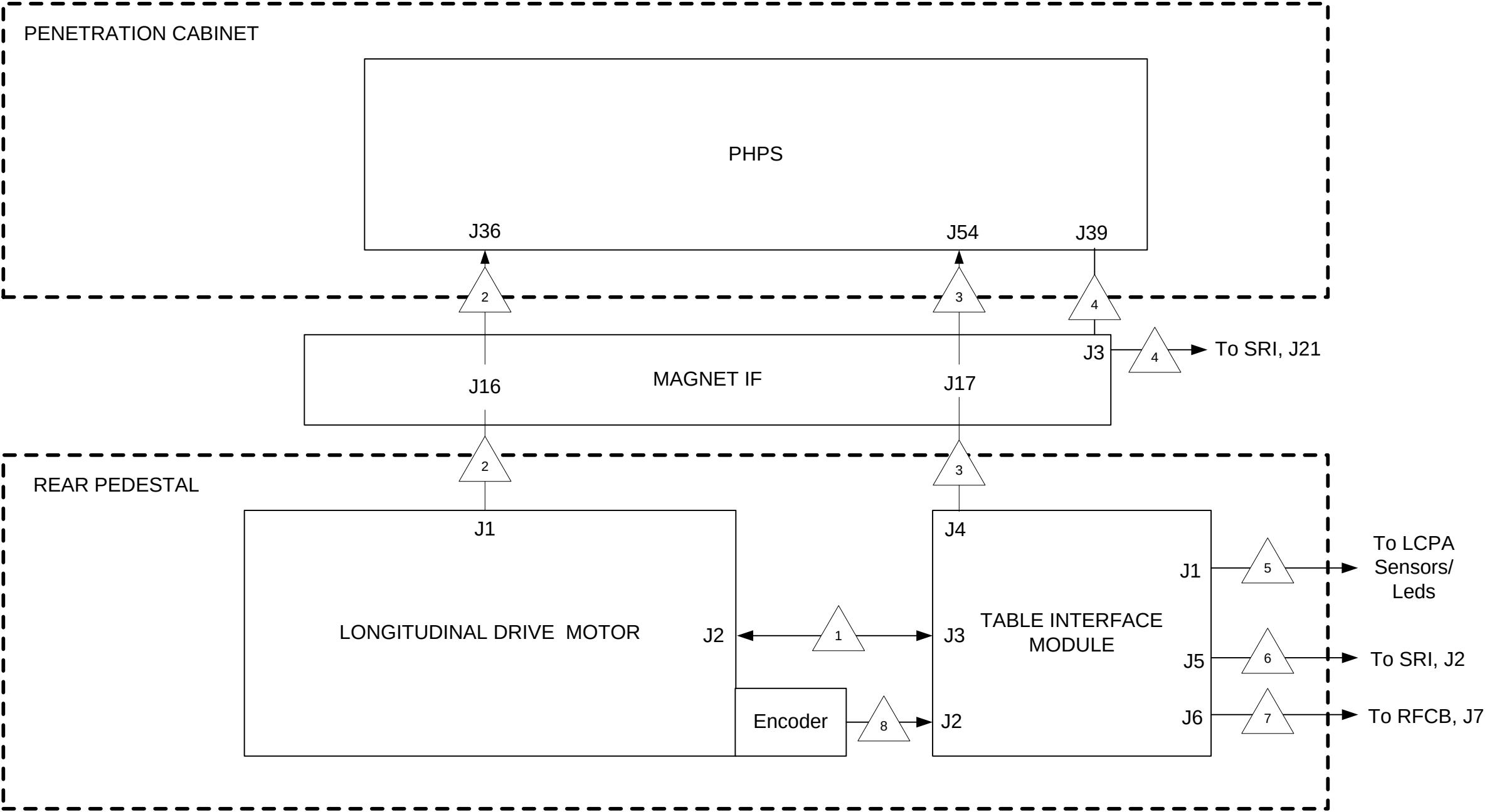


SIGNALS

- 1
- Control, Motor Encoder Channel A & B, Motor Encoder Index, Cable Detect for Motor Encoder, Reference for Temperature Signal, Digital 5V Power, Hall U, V, W Signal for Commutation, Motor Temperature Monitor
- 2
- Power, Motor Command to U,W, V Terminal, Long Drive Motor Cable Detect
- 3
- Hall U, V, W, Commutation Signals, Motor Temperature

- 4
- MCPDB Interface Connector, Dock Motor Sense, CAN RX Loopback to SRI4, Power Monitor Indicator, Motion Controller Amp Fault, Dock Power Cable Detect, Long Drive Cable Detect, Motor Encoders A & B, Motor Encoders A & B, Motor Encoder Index, Controls Dock Up Motor, Controls Dock Down, Fiber Optic Loopback Enable from SRI4, Motion Controller Amp Reset from SRI4, Motor Enable from SRI4, GND For MCPDB Cable Detect, Dock Motor Sense

- 5
- A Green, A Red, Home Limit Switch, Egress Limit Switch, 5V, EOT Limit Switch, LPCA Cable Detect
- 6
- Motor Encoder Signals A & B, Home Limit Switch, Egress Limit Switch, 8V, 1-Wire, EOT Limit Switch
- 7
- +5V DC, Green LED Port A, Red LED Port A, Auxiliary 1-8, Cable Connect, Ground
- 8
- Motor Encoder Signals A & B, 5V

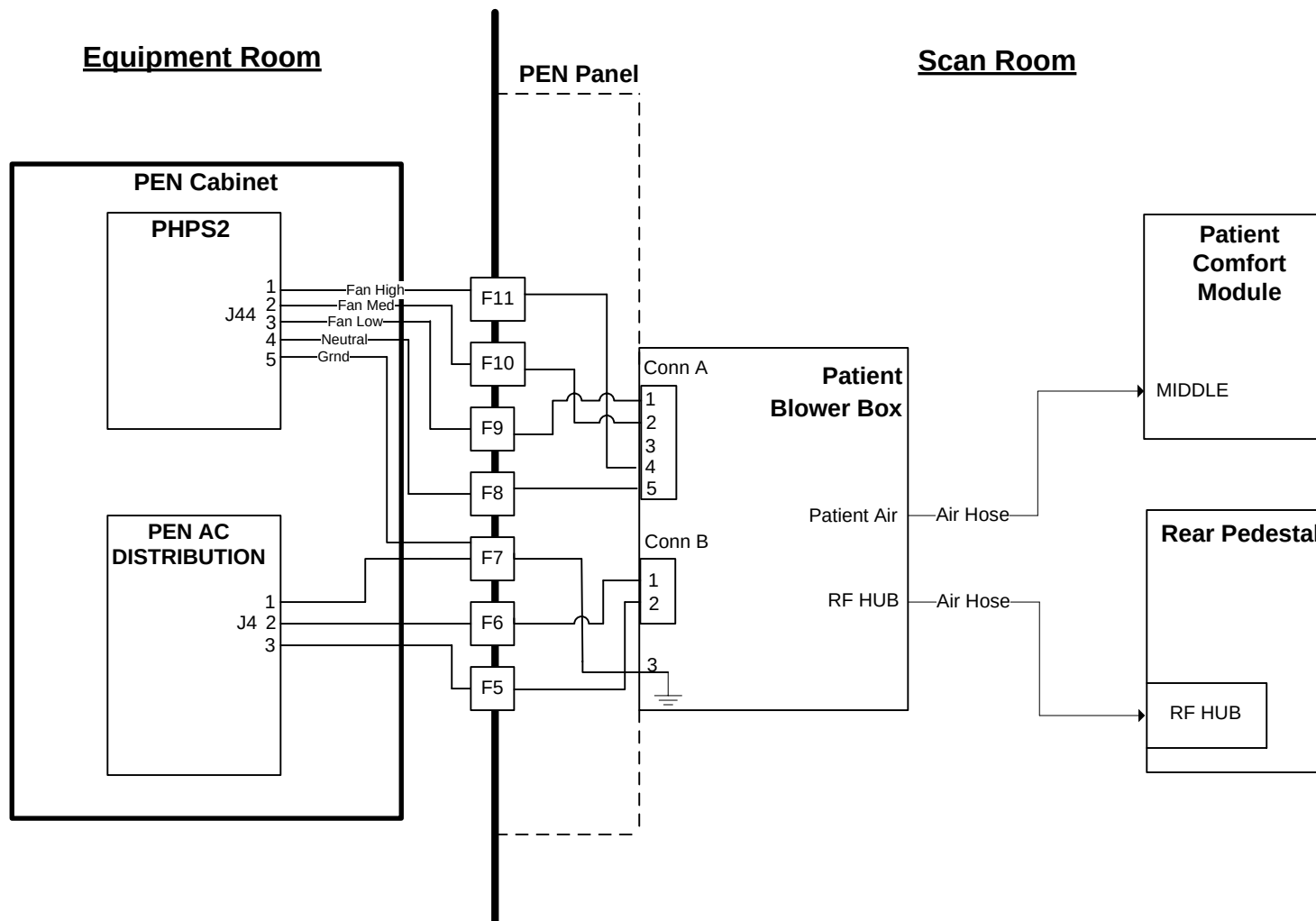**PATIENT HANDLING**

Patient Blower Box