



Power Supply Upgrade and Battery Replacement Installation Instructions GEHC 5376545

1. PROCEDURE

This procedure is applicable for replacing the 12Vdc battery and upgrading the EDS 11027 battery charger with a power supply charging circuit in the following panel:

GEHC P/N 5161287 Discovery MR750 Main Disconnect Panel which may also be referred to as the DVMR Main Disconnect Panel.

Reminder:



RISK OF ELECTRICAL HAZARD

Ensure all power sources feeding this MDP panel have been turned off and securely locked-out and tagged-out using the GEHC and hospital Lock-Out/Tag-Out (LOTO) procedures.

Depending on local rules, use protective equipment and tools. (Gloves and Arc Flash mask)

Required Tools & Equipment


- Screwdrivers, various sizes, slotted & Phillips
- Digital Voltmeter
- Wire cutter & wire stripper

Familiarize yourself with this entire Installation Procedure before proceeding.




2. POWER SUPPLY ASSEMBLY INSTALLATION

- a. Open the enclosure door.
- b. Ensure CB5 is off.
- c. Disconnect the wires connected to the battery.
- d. Locate the battery charger within the panel.
- e. Remove and discard the charger retaining bracket and its mounting hardware.
- f. Remove the charger from the duplex outlet.
- g. Cut the wires off of the charger as close as possible to the strain relief.
- h. Remove and discard the duplex outlet cover plate.
- i. Remove and discard the (2) screws that anchor the duplex outlet to the outlet box.

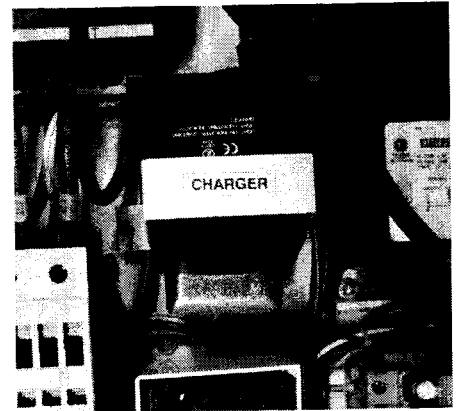


**CHARGER WIRES
(Cut from charger.)**

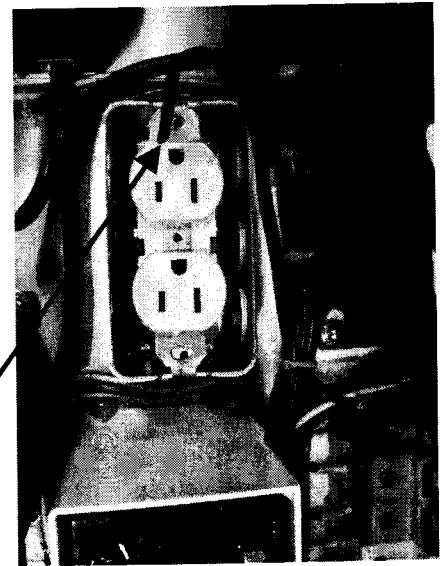

- a. Mount the Power Supply Assembly to the top of the outlet using the supplied #8-32 X 1/2" fillister head machine screws.

NOTE: It may be easier to mount the plate by removing the power supply from the DIN rail.


**Fillister head
machine screws**

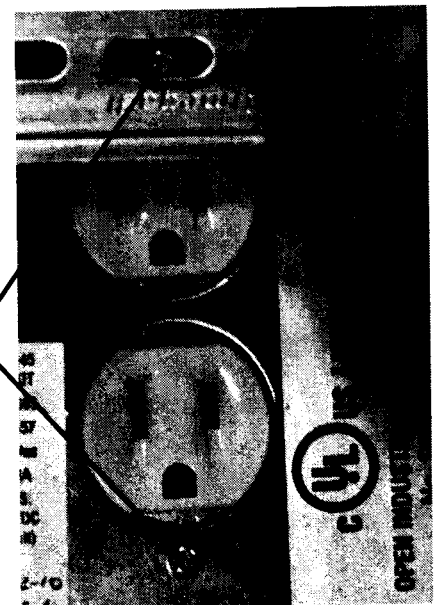




Install the power supply back on the DIN rail, adjusting the position to allow access other panel devices in the vicinity. Slide the end anchors tight against the power supply. Tighten the end anchors.

- j. Plug the power cord into the exposed receptacle.
- k. Open FU4 on the new power supply assembly.
- l. Locate the 2 conductor cord that was trimmed off of the old charger.
- m. Strip and connect the solid black conductor to the bottom side of terminal #12. Torque to 14.2 lb-in.
- n. Strip and connect the black w/ white stripe wire to the bottom side of the FU4 fuseholder. Torque to 20 lb-in.

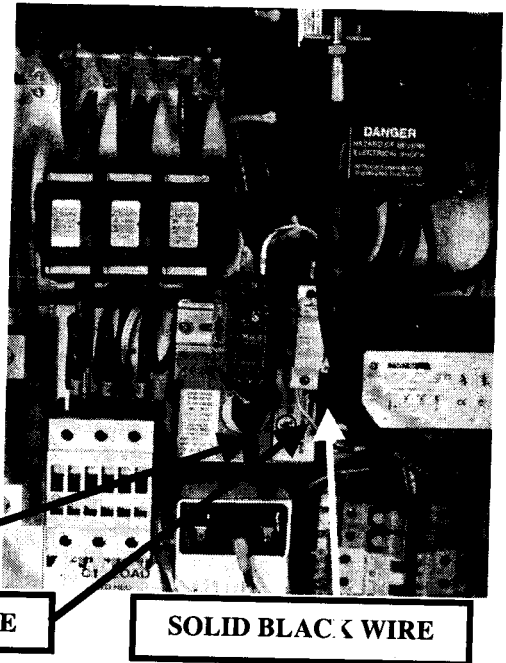

3. LABEL THE PANEL

- a. Apply the adhesive backed schematic drawing directly over the existing schematic drawing.
- b. Apply the PS1, FU4 label next to the schematic drawing.

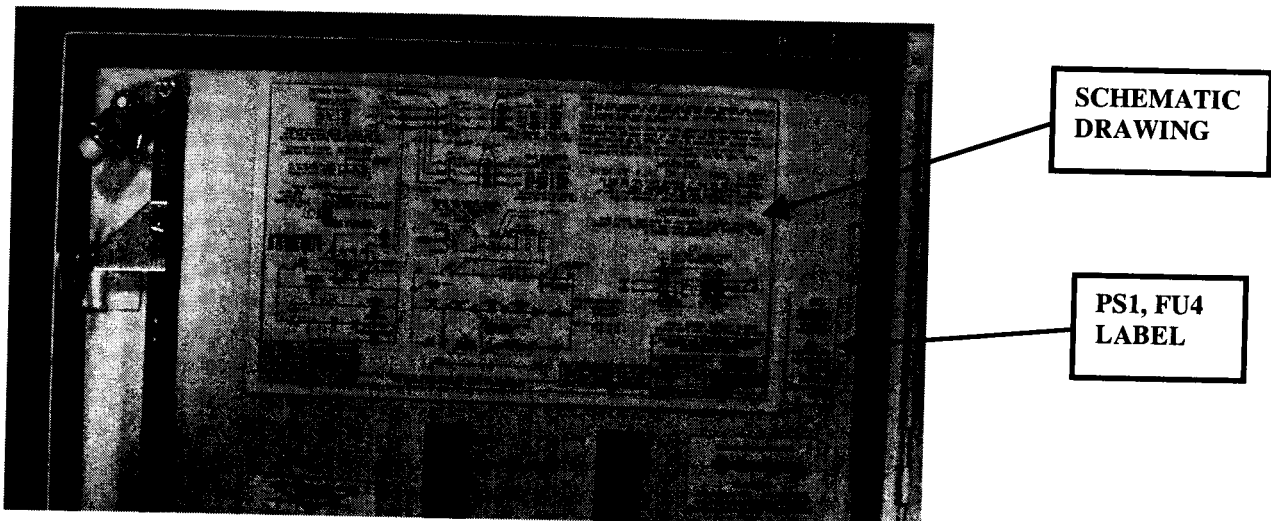

4. BATTERY INSTALLATION

- a. Refer to the installation instructions included with GEHC FRU 5390238.



5. VERIFY THE POWER SUPPLY OUTPUT VOLTAGE

- Ensure the FU4 fuseholder is open.
- While measuring the DC voltage between the load side of the power supply [wires #11A (+)] wire #12, (-), adjust the power supply output voltage for 13.45 to 13.55 VDC.
- Close the FU4 fuseholder.
- Monitor the battery terminal voltage for several minutes. Until the battery is fully charged, the terminal voltage will be increasing towards a final voltage of 13.45-13.55Vdc.

(The PULS P.S. green light should be on)

- Close the enclosure door and place the system back into operation.
- Dispose of the old power supply using proper GEHC disposal procedures.

6. COMPLETE AND RETURN THE ATTACHED FORM

Completing the attached form will allow Bevco Engineering to update our records confirming that this panel has been retrofitted with this upgrade kit.

Rev	Description of Change	Changed By	Approved By	Date
A	Initial Release	JFH		



DVMR Power Supply Upgrade Kit Installation

GEHC 5376545

Please complete this form and Fax it to Bevco Engineering.

Attention: Engineering Department

(262) 820-2599

Panel Serial Number: _____

Original Date of Mfg.: _____

Panel
Serial
Number

Original
Date of
Mfg.

Date of Upgrade Installation: _____

Upgrade Kit Serial Number: _____

Upgrade
Kit Serial
Number

