



GE Medical Systems

Technical Publications

Direction 2401983
Revision 4

Signa® IIQ kit Installation & Steering

Copyright © 2007 by General Electric Co.

Service Documentation

Jan 23, 2007

This document is intended to aid the Field Engineer in identifying the system type and required parts needed as part of the IIQ Kit installation process. Prior to installing the IIQ kit, the Field Engineer must determine magnet type and verify all the required kits and parts are on site.

Magnet Bridge Support Kit; P/N 5195386 / CAT M3090NE contains Bridge Supports for the following Magnets:


- S-series magnet enclosure parts are found in kits 5195382, 5195383 and 5195384
- LCC Wide Open (Illustration 1) enclosures parts are found in kits 5195383 and 5195384
- Cx / LCC Horizon (Salt Block) (Illustration 3) enclosure parts are found in kits 5195382, 5195383 and 5195384

The IIQ Kit is one of a number of kits that maybe required to properly upgrade the magnet to the HD / HDx current production level. When this IIQ upgrade is performed, it is highly recommended that the Field Engineer check the list and chart below for the scanner type and verify all the parts are on hand before beginning the upgrade.

- S-series Magnets
 - o RF Ground Strap – M3090NF
 - o Gradient Lead Board – 5198802 – ONLY when 2 AWG Leads BRM is used
- Cx Magnets
 - o RF Ground Strap M3090NF
 - o Gradient Isolation Kit - Coil Lift Tool 2261179 is required to install Gradient Isolation Kits
 - 2242380 BRM Gradient Isolation Kit
 - 2247493 CRM Gradient Isolation Kit
 - o Gradient Lead Board – 5220794 – ONLY when 2 AWG Leads BRM is used
 - o Gradient Cable Restraint Block 5182681-2 - ONLY when 2 AWG Leads BRM is used

Note : See the BRM/CRM upgrade chart below for details

BRM Upgrade Chart

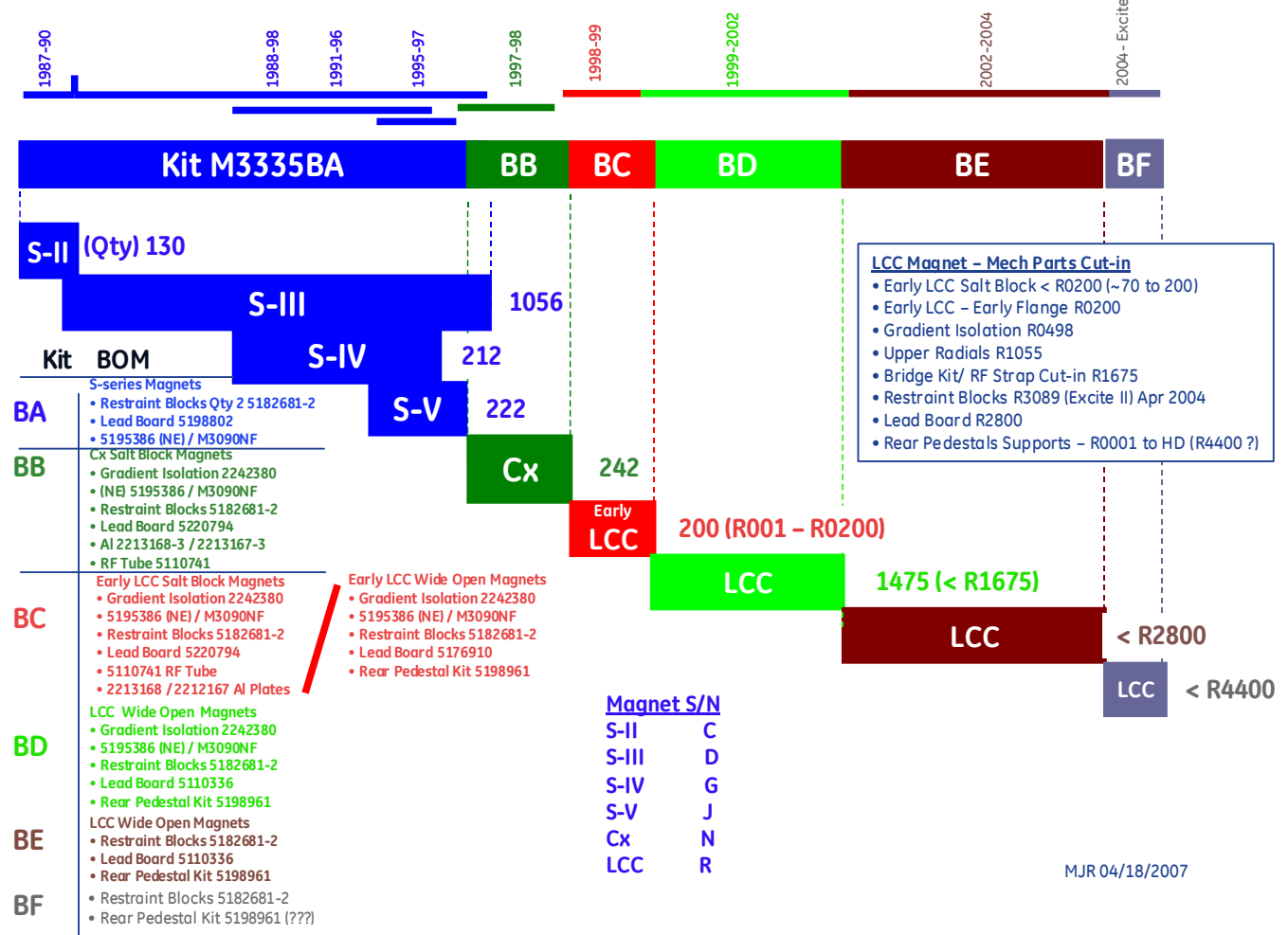

Wide Open Enclosure (Skintight)

- For Skintight (wide Open) covers, install kit labeled Bridge Support for LCC Magnet. Part # 5195383 and parts of 5195384 (Both included in IIQ kit 5195386). See Illustration 1. The installation document (2318581) is included in the kit. This installation document will provide further information on which parts are to be used from the kits.



Illustration 1

2. When the system has skintight covers, verify that Gradient Isolation brackets are installed. See Illustration 2. If the Gradient Isolation Brackets are missing, they must be installed.

- Kit 2242380 for BRM gradient coil
- Kit 2247493 for CRM gradient coil

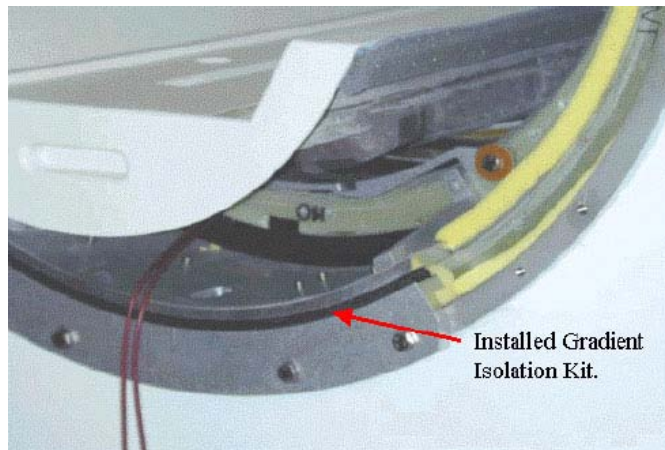


Illustration 2.

3. Important: Coil Lift kit (part number 2261179) must be ordered from the tool depot and on site before you attempt to install the gradient isolation brackets
4. Gradient Radial Supports kit (5195385) is included for magnets that need upper Radial support. (See direction 2318581 for more details).

Horizon (Salt Block)

1. For Horizon (Salt Block) covers, install kit labeled Bridge Support for S series (Long Bore), CX & LCC Magnet. Parts of kits 5195382, 5195383 and 5195384 (All included in IIQ kit 5195386). See Illustration 3. The installation document (2318581) is included in the kit. This installation document will provide further information on which parts are to be used from the kits.
2. Gradient Radial Supports kit (5195385) is included for magnets that need upper Radial support. (See direction 2318581 for more details)



Illustration 3