

Cradle Home Table Down Interlock Adjustment

([Click Here for Summary Version](#))

1 Personnel Requirements

Required Persons	Preliminary Reqs	Procedure	Finalization
1	Not Applicable	15 mins	Not Applicable

2 Overview

This procedure is applicable to the following table products:

- Signa Lightweight Patient Transport Table
- Signa Oncology Patient Transport Table Option
- Signa OR-Compatible Patient Transport Table Option

This adjustment is for the mechanical interlock which prevents the Patient Transport from being undocked or lowered while the cradle is in the magnet bore.

3 Procedure

1. Remove all upper and lower trim covers from sides of Patient Transport.
2. With cradle NOT HOME (removed), check the interlock mechanisms for tolerances shown in [Illustration 1](#). To adjust:
 - a. 1/16 to 1/32 inch (1.59 to 0.79 mm) gap - carefully bend interlock rod. See [Illustration 1](#).
 - b. 1/8 inch (3.18 mm) gap - loosen lock nut and tighten or loosen the cable adjustment screw to achieve the proper gap. Tighten lock nut. See [Illustration 2](#).

Illustration 1: CRADLE HOME/TABLE DOWN INTERLOCK

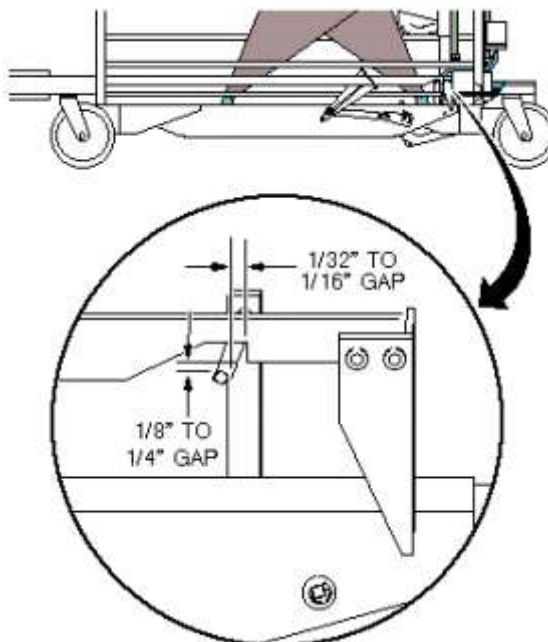
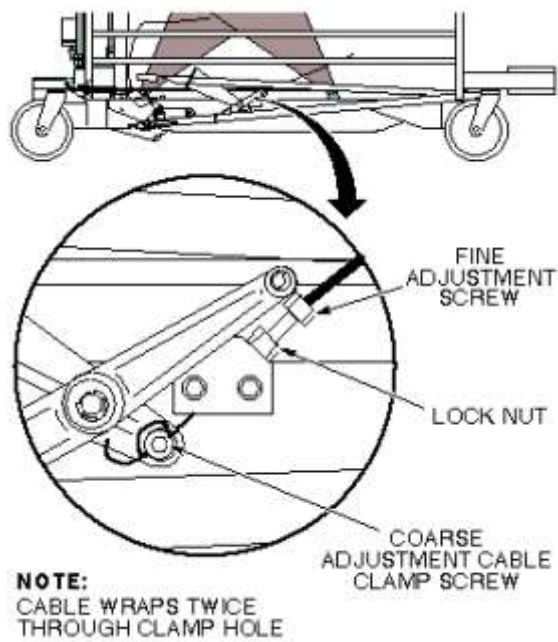


Illustration 2: CRADLE HOME/TABLE DOWN INTERLOCK

NOTE: If the adjustment screw does not provide enough range for adjustment, coarse adjustment may be made to the cable length where the cable end is clamped into actuator bracket. Fine adjustment will be required after performing a coarse adjustment.

4 Finalization

No finalization steps.

[Back to top](#)