

GE Healthcare MR Class A eService Notes are copyright (c) to General Electric Company, Inc. with all rights reserved.
Duplication and distribution of any class of MR eService Notes, in any form, by entities other than GE Healthcare is strictly prohibited.

SERVICE NOTES

Feb 15, 2012

Class A

Discovery MR750w systems

SN DOC1102422 Rev1

Discovery MR750w Service Methods: Error Correction of Body Coil Tuning Procedure

EFFECTIVITY: Discovery MR750w Service Methods
(Class M: 5670006 Rev1, Class A: 5670008–3 Rev1)

PROBLEM: Two errors were found in Body Coil Tuning of Discovery MR750w Service Methods Rev1.
(Adjustments/Calibrations / Radio Frequency (RF) / Body Coil Tuning)
These errors were reported by SPR (HCSDM00117041 and HCSDM00119031).

1. Phantom Setup in 4.2 Hardware Setup / Step3 is wrong.
2. J numbers of TRSW in "Illustration 14: Hardware Setup" is wrong.


SOLUTION:

Solution for Problem#1:

The front edge of body loader should be set to **835mm (33")** away from the front edge of the Bridge to align it to the Magnet center. However, the procedure shows wrong value "990mm (39")" in Service Methods Rev1. The description in "Body Coil Tuning / 4.2 Hardware Setup / Step3" should read as follows.

3. Place Body Loader with phantom on to the left-right center of the bridge and push it to the Magnet Center (The front edge of body loader is **835mm (33")** away from the front edge of the Bridge).

Also, Illustration 13 should be corrected as follows.

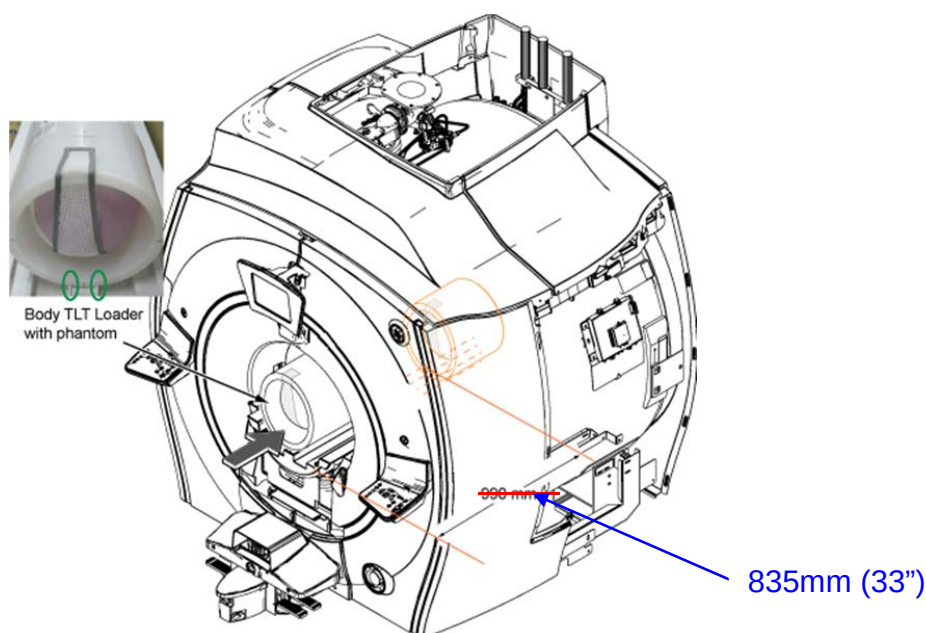


Copyright© 2012 by General Electric Company

GE Healthcare MR Class A eService Notes are copyright (c) to General Electric Company, Inc. with all rights reserved.
Duplication and distribution of any class of MR eService Notes, in any form, by entities other than GE Healthcare is strictly prohibited.

Solution for Problem#2:

J numbers of TRSW in Illustration 14 should be corrected as follows.

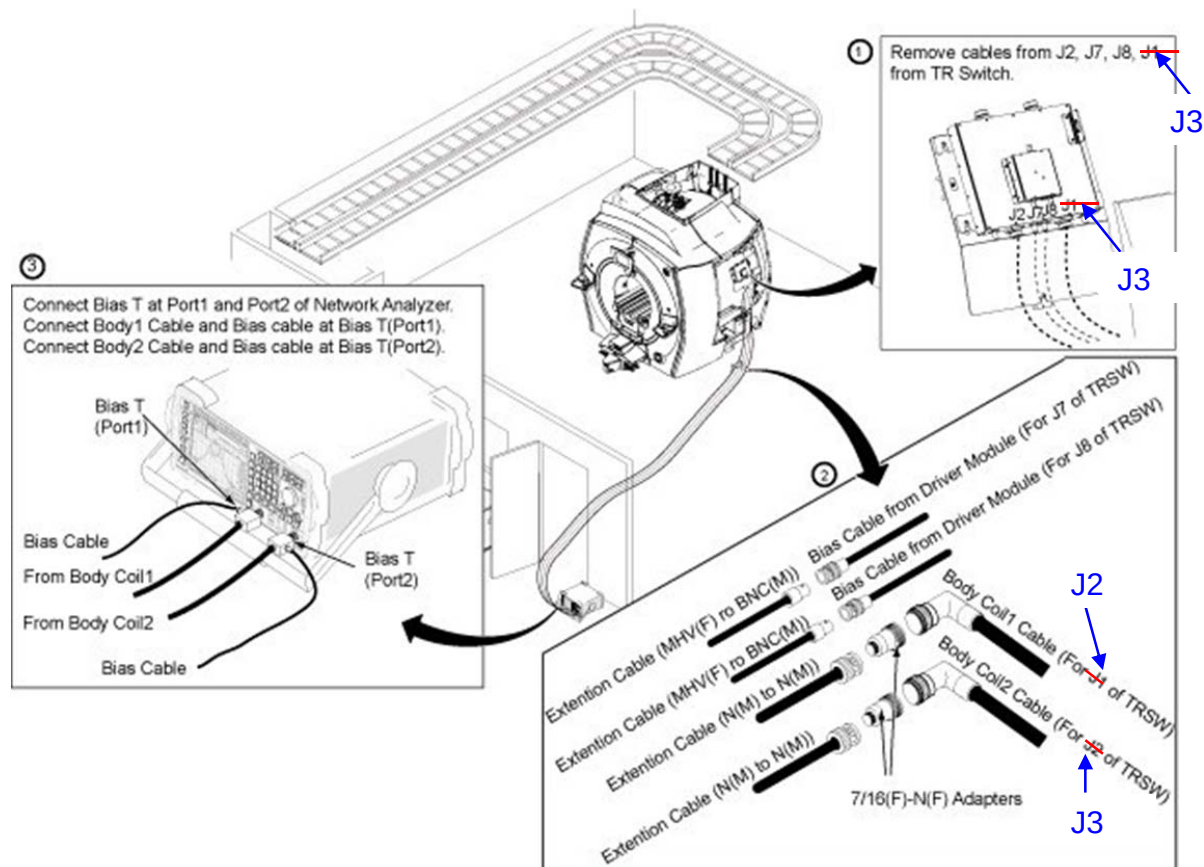

These errors will be corrected in the next revision of Discovery MR750w Service Methods.

Next revision: Class M: 5670006 Rev2, Class A: 5670008-3 Rev2

AUTHOR: Yutaka Masumo
GE Healthcare Japan - MRIB/SVD Gr.