

MAGNETOM *syngo* Harmony/Symphony/Sonata, Concerto, Trio, Allegra

MR

Installation Instruction

Installation of external MOD for XP-based Systems

Revision: 01

Date: 28.10.2004

Name:
Dept.:

Schülein
CS PS 23

© Siemens AG 2004
The reproduction, transmission or use
of this document or its contents is not
permitted without express written
authority. Offenders will be liable for
damages. All rights, including rights
created by patent grant or registration
of a utility model or design, are
reserved.

Doc. Gen. Date: 10.04

General

Describes how to install an external MOD drive on Windows XP-based MR-systems

Prerequisites


System is running NUMARIS/4 VA25A or 4VB11A/D

Installation

For systems based on Windows XP R610 PCs are used which are equipped with a SCSI-controller. Below the backside of the PC with all the connections are described.

- Remove the SCSI - terminator (4)
- Connect the SCSI cable of the external MOD drive instead

MRC:

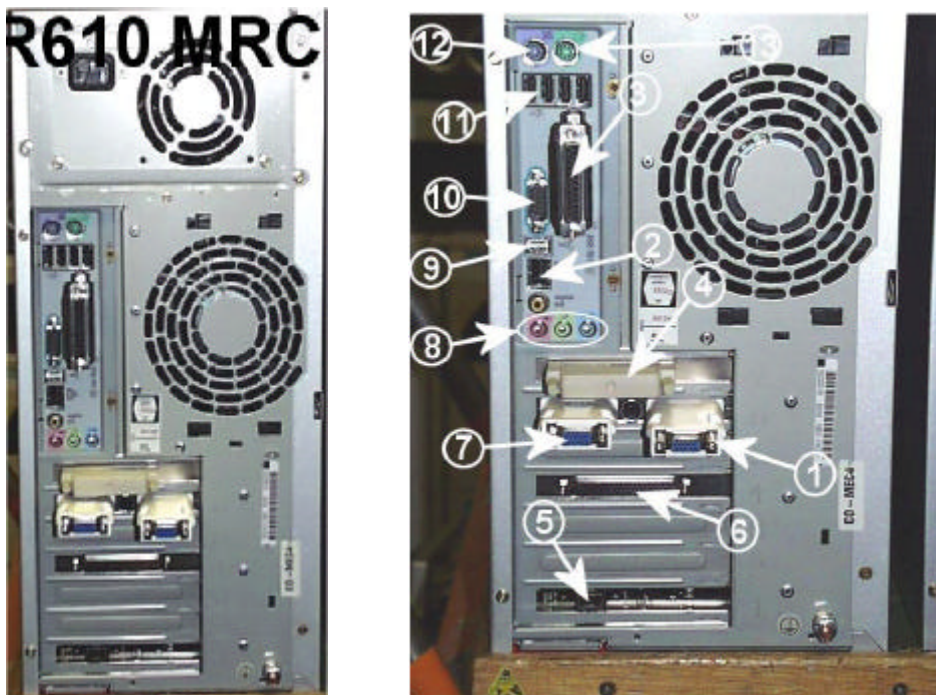

1. For the LCD display (VGA Output)
2. For the measurement network (LAN connector onboard). Valid only for software version syngoMR2002B.
3. For the dongle and printer (parallel connection)
4. For external SCSI devices like the MOD drive (remove the SCSI terminator and connect the external MOD drive if this option is available)
5. For the hospital network (additional LAN card). Valid only for software version syngoMR2002B.
6. For the serial port extender (MOXA card, system standby function)
7. Additional VGA output for a second monitor
8. Line In Jack (blue)- Connects to external devices for playback or recording
Headphone Out Jack (green) - Connects to headphone or speakers (Amplifier integrated)
Microphone In Jack (pink) - Connects to an external microphone
9. Firewire connector
10. Serial connector
11. USB 2.0/1.1 connectors
12. Keyboard connector
13. Mouse connector

- Connect the power cable of the MOD drive
- Start up the system

Configuration of the external MOD Drive

Perform the following steps logged in as **meduser**:

- 1) In the **Options** pull-down menu, select the **Local Service**.
- 2) Enter the valid service key -> and click **OK**.
- 3) Select **Configuration**.
- 4) Accept the Options page by selecting **Next**.
- 5) Select **offline devices**
- 6) Enter MOD for 'Logical Name'
- 7) Select the radio button 'Optical Disk Drive' -> the drive letter will be determined by the program (Usually F: if nothing else is connected to the system.)
- 8) Check the **SPI file system** box for the SCSI settings.
- 9) Click **Save**.
- 10) Click **Finish**.
- 11) Reboot the computer.

Final test

Check if the MOD appears in the patient browser.

This page intentionally left blank.