

SECTION 6 - MOVING MAGNET TO MR SUITE

WARNING!

LOSS OF MAGNET VACUUM WILL RESULT IN THE VIOLENT EXPULSION OF HELIUM GAS, WHICH CAN CAUSE ASPHYXIATION IN ENCLOSED AREAS. USE EXTREME CAUTION TO NOT CONTACT / DAMAGE THE VACUUM VESSEL DURING MAGNET TRANSIT OR SITING.

IMPORTANT !!!

Make sure that Magnet Mounting Plate has been installed in Magnet Room and that magnet orientation is determined prior to moving magnet.

WARNING!

DO NOT APPLY ANY FORCE TO FIBERGLASS COVERS ON TOP, SIDES AND ENDS OF MAGNET.

DO NOT APPLY FORCE TO THE SINGLE-STAGE COLDHEAD SHOWN IN ILLUSTRATIONS 6-1, 6-2 AND 6-3.

USE LIFTING STRAPS / CABLES AS SHOWN TO LIFT MAGNET (REFER TO SECTION 3-5).

Note

Identify the magnet location in the Magnet Room and check all clearances along the magnet delivery route before magnet delivery. See Illustrations 6-1 and 6-2.

1. Use the Lifting Rails pre-installed on the magnet for jacking to raise the magnet load points for roller dolly installation. Use shims / plates when rolling magnet with roller dolly over door thresholds and other inclines. See Illustrations 6-3.

Note

Roller dollies are recommended for moving magnet into the Magnet Room as shown in Illustration 6-3. If roller dollies are used, place steel floor plates along the magnet delivery route, as well as over door thresholds and other inclines.

SECTION 6 - MOVING MAGNET TO MR SUITE (continued)


2. Move magnet to Magnet Room. If using a motorized tow vehicle, attach straps / cables / chains to Lifting Rails with shackles. Do not apply any loads to any enclosure cover parts nor allow straps / cables / chains to scrape enclosure cover parts.

Note

If there are turns in the delivery route, adjust roller dollies to appropriate positions to negotiate turns. The height can be reduced by putting roller dollies under Lifting Rails.

3. If there is not enough vertical clearance for the magnet to clear headers or doorways, call the MAC Team representative. The height can also be reduced by putting roller dollies under the lifting rails.
4. Identify table end of magnet before moving into Magnet Room. Continue with Section 7, Magnet Mounting.

SECTION 6 - MOVING MAGNET TO MR SUITE (continued)

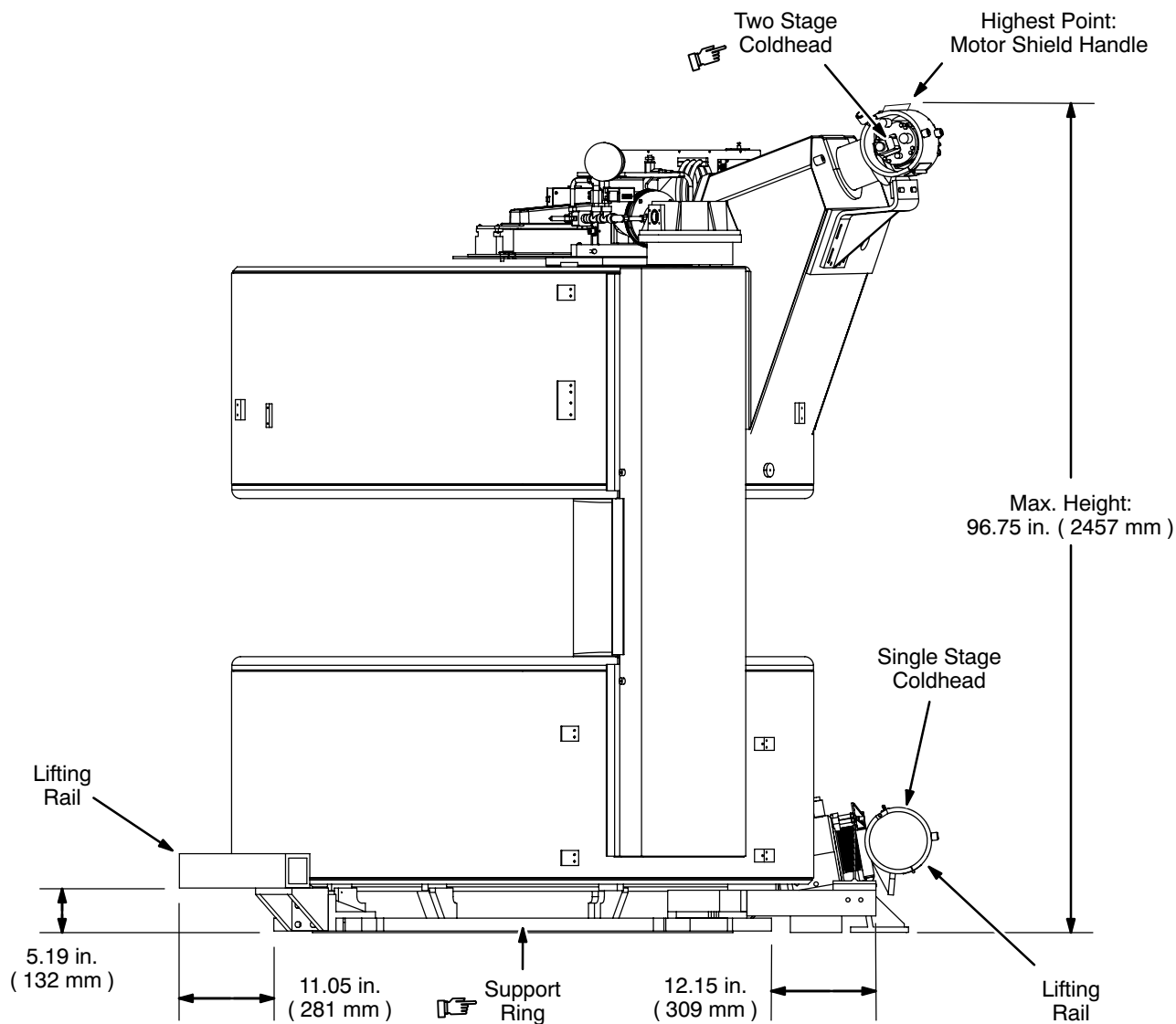


MAGNET CLEARANCES
ILLUSTRATION 6-1

Note

Minimum height for service clearance is 102.68 inches (2608.1 mm). This clearance is needed for ramp lead, shim lead and fill line installation.

SECTION 6 - MOVING MAGNET TO MR SUITE (continued)

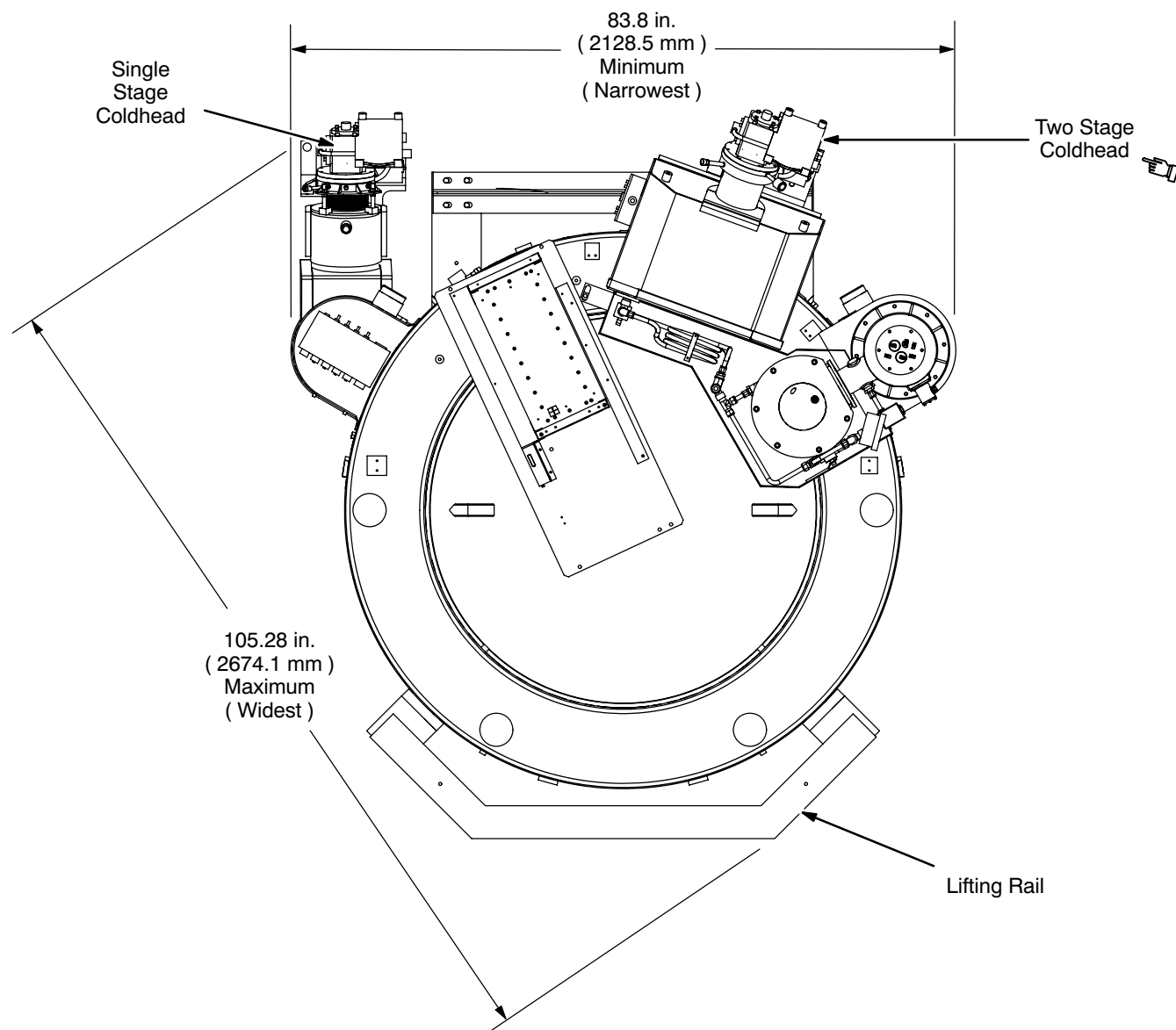

MINIMUM / MAXIMUM MAGNET DIMENSIONS
ILLUSTRATION 6-2

SECTION 6 - MOVING MAGNET TO MR SUITE (continued)**CAUTION**

Keep magnet level at all times. Uneven jacking of corners could result in magnet shifting on jacks.

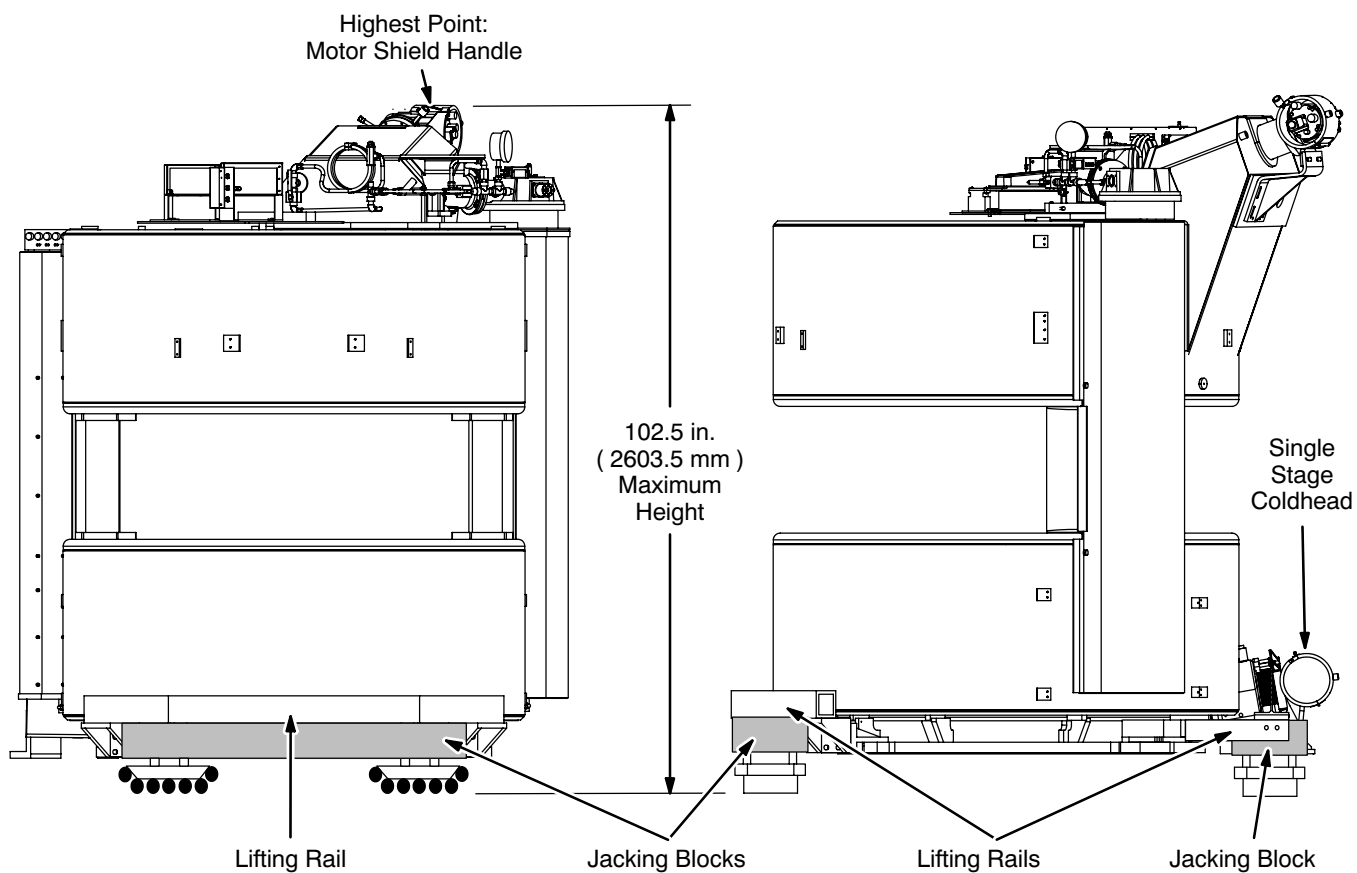
**MAGNET WITH ROLLER DOLLIES ATTACHED**

ILLUSTRATION 6-3

SECTION 6 - MOVING MAGNET TO MR SUITE (continued)