





1.5T SIGNA MR/i – W/EXCITE HD

RECOMMENDED MAGNET ROOM
CEILING HEIGHT 9'-1" [2768mm]

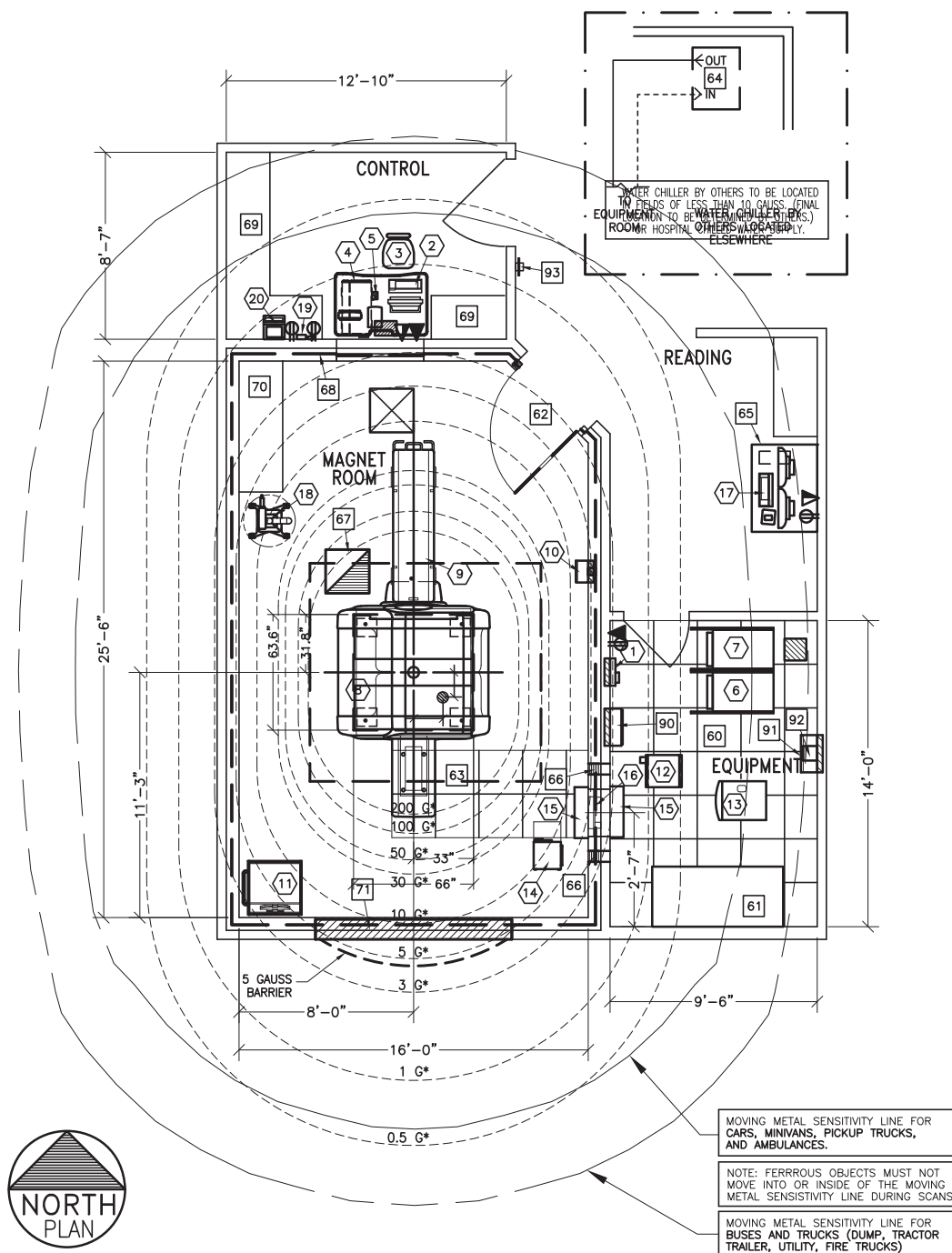

REFER TO PRE-INSTALLATION MANUAL FOR ROOM SIZES BASED UPON CONFIGURATION

ROOM DIMENSIONS	LENGTH & WIDTH	CEILING HEIGHT
MIN. MAGNET RM:	19'-8" X 11'-1" [5.99m X 3.39m]	8'-9" [2.67m]
MIN. EQUIP. RM:	8'-0" X 12'-0" [2.44m X 3.67m]	8'-0" [2.44m]

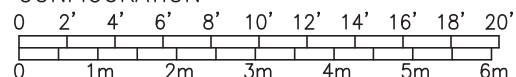

SEE APPENDIX "A" FOR POWER SPECIFICATIONS

ROOM PLAN LOCATED ON PAGE 8-185

ITEM	QUANTITY	ITEM DESCRIPTION	WEIGHT	HEAT OUTPUT
①	1	MAGNET MONITOR	22 lbs	204 btu
②	1	OPERATOR WORKSPACE W/COLOR LCD MONITOR	125 lbs	4948 btu
③	1	OPERATOR'S CHAIR		
④	1	OPERATOR WORKSPACE CABINET	198 lbs	
⑤	1	PATIENT ALERT CONTROL BOX		
⑥	1	RFS CABINET	780 lbs	27071 btu
⑦	1	HFD/PDU CABINET	1805 lbs	34129 btu
⑧	1	1.5 TESLA LCC ACTIVE SHIELD MAGNET	12608 lbs	8191 btu
⑨	1	PATIENT TRANSPORT TABLE (DOES NOT INCLUDE PATIENT)	279 lbs	
⑩	1	MAGNET RUNDOWN UNIT	8 lbs	
⑪	1	SPT PHANTOM CABINET	350 lbs	
⑫	1	SHIELD COOLER CABINET	275 lbs	170 btu
⑬	1	WATER CHILLER FOR BRM BODY COIL	295 lbs	13993 btu
⑭	1	BLOWER BOX	46 lbs	3412 btu
⑮	2	PENETRATION PANEL COVERS		
⑯	1	RF PENETRATION PANEL	88 lbs	324 btu
⑰	0	ADVANTAGE WORKSTATION WITH TWO LCD MONITORS (OPTIONAL)	81 lbs	1109 btu
⑱	0	INJECTOR HEAD ON PEDESTAL (OPTIONAL)	59 lbs	
⑲	0	BATTERY CHARGING UNIT (OPTIONAL)	4 lbs	
⑳	0	CONTROL ROOM UNIT (OPTIONAL)	15 lbs	

ENVIRONMENT

Ambient Operation Temperature: 50-70° F [10-21° C]
 Allowable Temperature Change: 5° F [3° C] per hour
 Humidity: 30-60 per cent non condensing
 Allowable Humidity Change: 10 per cent per hour
 Max - Altitude: 8000 ft. [2438m] above sea level

ROOM PLAN LOCATED ON PAGE 8-185

ITEM 	QUANTITY	ITEM DESCRIPTION	WEIGHT	HEAT OUTPUT
60		***** CUSTOMER/CONTRACTOR SUPPLIED AND INSTALLED ITEMS *****		
61		ACCESS FLOOR WITH 2' x 2' (610 x 610mm) REMOVABLE PANELS.		
62		AIR CONDITIONING. (VIBRATION ISOLATION IS RECOMMENDED AT SUPPORTS OF EACH UNIT TO BE INSTALLED.)		
63		MINIMUM DOOR OPENING FOR EQUIPMENT DELIVERY IS 43 IN. W x 82 IN. H [1092mm x 2083mm], CONTINGENT ON A 96 IN. [2438mm] CORRIDOR WIDTH		
64		NON-METAL ACCESS FLOOR WITH 2' x 2' (610 x 610mm) REMOVABLE PANELS & SUPPORT HARDWARE REQUIRED WITHIN MAGNET ROOM.		
65		WATER CHILLER		
66		WORKSTATION TABLE		
67		RF FILTERS - LOCATE WITHIN 24 in. [610 mm] OF THE RF COMMON GROUND STUD		
68		MAGNET ROOM EXHAUST FAN		
69		RF SCREEN, INCLUSIVE OF WALLS, FLOOR, DOOR, ETC. GROUND IMPEDANCE GREATER THAN 1000 OHMS. ATTENUATION 100dB AT 10-100MHZ PLANEWAVE.		
70		COUNTERTOP WITH DRAWERS FOR MISCELLANEOUS ITEMS.		
71		BASE CABINET FOR STORAGE OF: SURFACE COILS, PATIENT POSITIONING PADS, PHANTOMS, ETC.		
90		MINIMUM 9 FT. -0 IN. [2743 mm] x 9 FT. -0 IN. [2743 mm] REMOVABLE WALL SECTION FOR MAGNET DELIVERY/REMOVAL.		
91		MAIN DISCONNECT CONTROL - 94 lbs. (43 kg.), 900 BTU/HR (264W) - CAT NO. E4502SP FOR 480-3 WYE.		
92		DC LIGHTING CONTROL PANEL 155 lbs (70 kg) 1024 BTU/HR. (CAT. NO. E4502SC/SE - BASIC SYSTEM)		
93		DC LIGHTING AUTO TRANSFORMER 60 lbs [27 kg] (PART OF VARIABLE DIMMER SYSTEM) (CAT. NO. E4502SD/SF INCLUDES BASIC SYSTEM)		
93		METAL DETECTOR (HAND HELD)		