

SECTION 5 – MOVING MAGNET TO EXAM ROOM

Description

Before moving magnet to exam room (fixed sites), make sure mounting holes, anchors and washers are installed as per Magnet Bolt Down Kit instructions. For normal installations use Magnet Bolt Down Kit P/N 2104219. For Seismic mounting use Magnet Bolt Down Kit P/N 2105959.

Note

Roller Dollies are recommended for moving magnet into the examination room as shown in Illustration 5–1. If Roller Dollies are used, place steel floor plates along the magnet delivery route.



Use the Jacking Pads located on the magnet to raise the magnet for Roller Dolly installation.

Use shims when rolling magnet with Roller Dolly over door thresholds and other inclines.

Note

The magnet location in the examination room must be marked. Use tape to mark the four corners of magnet on the examination room floor. Refer to architectural drawings to determine the exact location of magnet within the examination room.

1. Push magnet to exam room. If using a motorized tow vehicle, attach chains around the magnet base support pads and pull magnet.
2. If there are turns in the delivery route, adjust Roller Dollies to appropriate positions to negotiate turns.
3. Move Magnet into the Exam Room. If there is not enough vertical clearance for the Magnet, remove the 1.5 inch plates from the bottom of the magnet or the Vent Adapter from the Vertical Penetration to obtain additional clearance. If more clearance is still needed call the MAC Team Representative for further removal of Vertical Penetration components.
4. Position magnet on the examination room floor. If ceiling is in place, drop a plumb bob from the center of the vent hole and mark the floor at the location where the plumb bob points. This corresponds to the location of the center of the Vent Adapter on the Magnet.

Note


Table end of magnet (See Introduction, Illustration 4–1) must be oriented away from rear wall of exam room.



Keep magnet level at all times. Uneven jacking of corners could result in shifting of magnet on jacks.

5. Jack the magnet up sufficiently, at the 4 jacking pads, and remove the 4 Roller Dollies.
6. Slowly lower magnet onto the exam room floor.


ALL DIMENSIONS ARE IN INCHES (MILLIMETERS)
* DIMENSIONS INCLUDE 1.50 (38.1) PLATES ON BOTTOM OF EACH FOOT.

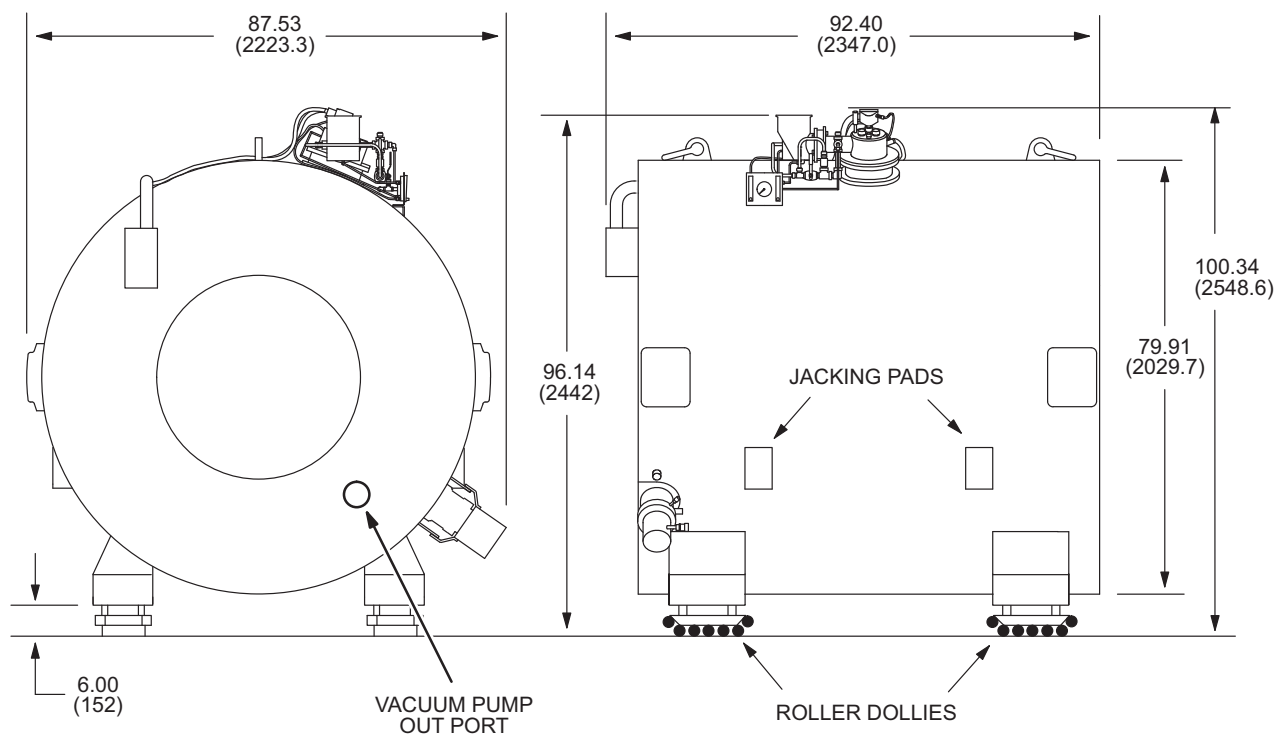

MAGNET WITH ROLLER DOLLIES ATTACHED
ILLUSTRATION 5-1

ALL DIMENSIONS ARE IN INCHES (MILLIMETERS)

* DIMENSIONS SHOWN WITHOUT FOOT BLOCKS ATTACHED TO MAGNET FEET.

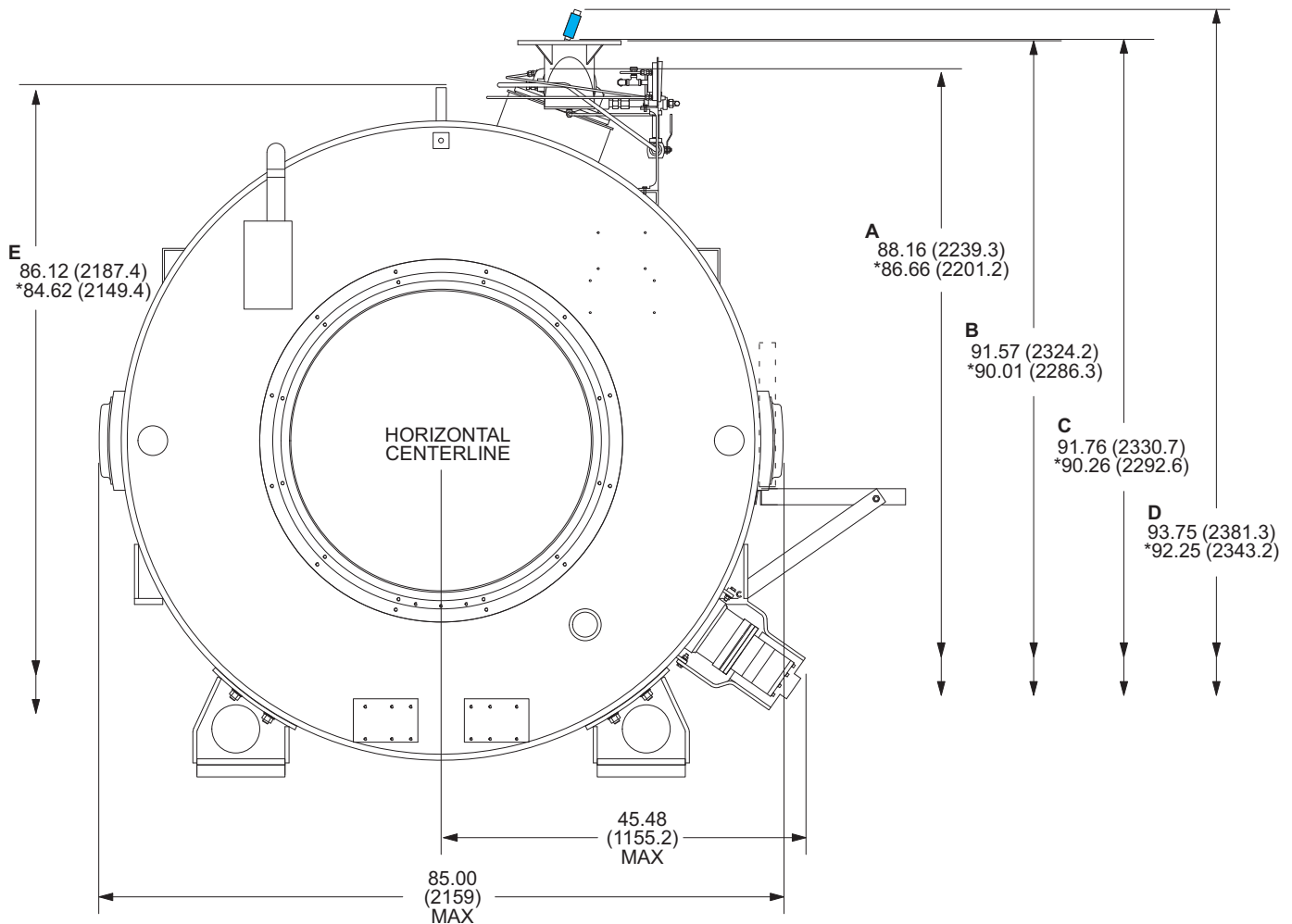

- A: DISTANCE TO UPPERMOST POINT ON TURRET COVER PLATE.
 B: DISTANCE TO TOP OF VENT ADAPTER.
 C: DISTANCE TO TOP OF SHIM LEAD CONNECTOR BOX.
 D: DISTANCE TO TOP OF SHIM LEAD BRASS KNURLED NUT.
 E: DISTANCE TO TOP OF LIFTING LUG

Note

Minimum height for service clearance is 102.68 inches (2608mm). This clearance is needed for Ramp Lead, Shim Lead, and Fill Line installation.

MAGNET CLEARANCE DIMENSIONS

ILLUSTRATION 5-2