

GE HealthCare
Technical Bulletin



Technical Bulletin Unique Identifier	DOC2943830
Title	Service Note of SIGNA MR355, SIGNA MR360, Brivo MR355 and Optima MR360 Magnet Enclosure Localization
Subject	This Service Note introduce the Magnet Enclosure of SIGNA MR355, SIGNA MR360, Brivo MR355 and Optima MR360 do Material/Manufacturing Method change from PVC to PC/ABS, with beam localization, for cost reduction purpose.
Applicability (e.g., affected product, process)	1.5T SIGNA MR355 1.5T SIGNA MR360 1.5T Brivo MR355 1.5T Optima MR360
Date	15-March-2024
Relevant Product Business (Select one)	☒ Imaging ☐ Patient Care Solutions ☐ Ultrasound ☐ Other: _____
Originating Modality	MR
Details	This Service Note introduce the Magnet Enclosure of SIGNA MR355, SIGNA MR360, Brivo MR355 and Optima MR360 do Material/Manufacturing Method change from PVC to PC/ABS, with beam localization, For cost reduction purpose.
Tracking Number	ECR2372347
Affected Service Publications	• SIGNA MR355/MR360 Service Methods (M:5856347-3 EN Rev: 7, C:5856355-3EN EN Rev: 8, A:5856356-3EN EN Rev: 7) • Optima MR360/Brivo MR355/Optima MR360 Advance/Brivo MR355 Inspire Service Methods (SV23.0 platform & SV23.1 platform) (M:6336165-3EN Rev: 4, C:6336163-3EN Rev: 4, A:6336161-3EN Rev: 4)
Resolution	See attached
Technical Bulletin Owner	SONG Yang (SSO: 503386740)

(This SN is classified as Class A)

GE HealthCare

Technical Bulletin



1. Effectivity:

1.5T SIGNA MR355

1.5T SIGNA MR360

1.5T Brivo MR355

1.5T Optima MR360

2. Overview:

This service note introduces that the magnet enclosures and beam of SIGNA MR355/SIGNA MR360/Brivo MR355 and Optima MR360 have changed supplier for cost optimization.

The enclosure Material/ Manufacturing Method change from PVC to PC/ABS.

The overall size, hardness, weight, outlook, mounting points, number of parts and installation process of enclosures and beams remain the same.

For FRU: legacy enclosure will be installed for the customer until the legacy enclosure is out-of-stock, after which the new enclosure will be installed.


For forward production: the package change from G50002AA to G50022CH+G50022DJ.

Scope:

Enclosures as shown in Figure 1 with Table 1 below use the new material (can mix use with legacy enclosure).

And “Beam Collector 450w GEM” part number change from 5415854 to 5415854-100 for localization purposes. As shown in Table 2. The overall size, hardness, weight, outlook, mounting points, number of parts and installation process remain the same.

Figure 1: The enclosures have changed.

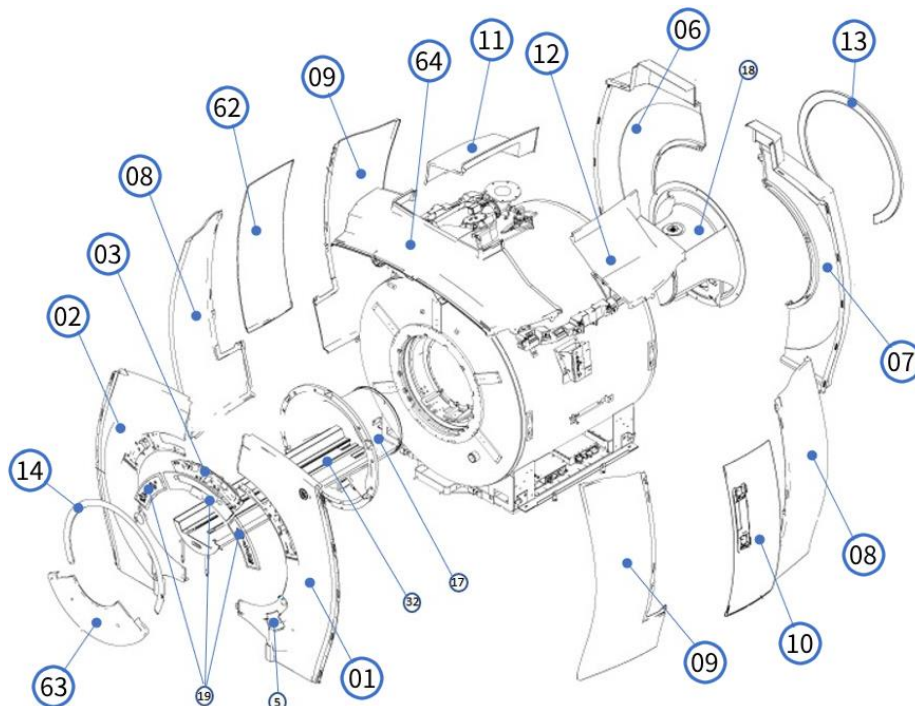

GE HealthCare

Technical Bulletin



Table 1: Enclosures use the new material.

NO.	Legacy Part No.	New Part No.	Description	Qty	FRU?	Gcat#
1	7500100	7500100-100	Cover Front RH Mulan Dual Port Locations-Asahi	1	Yes	G50022CH
2	5447195	5447195-100	COVER FRONT LH ASM 16 BEAT-Asahi	1	Yes	
3	5444782-2	5444782-102	COVER FRONT TOP ASM 16BEAT-Asahi	1	Yes	
6	5450848	5450848-100	Cover Rear LH Assembly 16 BEAT-Asahi	1	Yes	
7	5441557	5441557-100	Cover Rear RH Assembly 16 BEAT-Asahi	1	Yes	
11	5404661-2	5404661-102	COVER TOP, REAR, LH, BONDED, 16 BEAT-Asahi	1	Yes	
12	5404663-2	5404663-202	COVER TOP, REAR, RH, BONDED, 16 BEAT-Asahi	1	Yes	
13	5150411-2	5150411-102	Rear Trim Ring Assy for 16 BEAT-Asahi	1	Yes	
63	5161159	5161159-100	Front Bridge Cover Assembly-Asahi	1	Yes	G50022DJ
8	5404656-2	5404656-110	COVER SIDE LH BONDED,Signa High-Asahi. For 16 beat	2	Yes	
9	5404657-2	5404657-110	COVER SIDE RH BONDED,Signa High-Asahi. For 16 beat	2	Yes	
10	5457457	5457457-110	COVER SIDE CENTER RH ASM,Signa High-Asahi. For 16 Beat	1	Yes	
14	5150410-2	5150410-102	Front Trim Ring Assy for 16 BEAT-Asahi	1	Yes	
62	5404655	5404655-110	Side Cover, LH CTR, Bonded, Signa High-Asahi. For 16 Beat	1	Yes	
64	5404658	5404658-100	COVER, TOP FRONT ASM 450w-Asahi	1	Yes	

Table 2: Beam Collector 450w GEM.

Legacy Part No	New Part No	Description	FRU?	Gcat#
5415854	5415854-100	Beam Collector 450w GEM-Asahi	No	G50022CH

3. Pre- installation PIM

Enclosure with package shipping specifications:

Enclosure dimension still the same. And the package change from one kit G50002AA to two kits G50022CH+G50022DJ.

Both two new kits are smaller in size than the original kit, and their overall weight is reduced.

Gcat#	Type	Length x width x height (mm)	Weight (Kg)
G50002AA	Legacy	2680x1220x1510	316
G50022CH	New	2440x1540x1060	182
G50022DJ	New	2240x1020x1030	124

4. Installation

The overall size, hardness, weight, outlook, mounting points, number of parts and installation process remain the same. No service impacts.

5. Spare parts (FRU)

Replacement

The replacement process remains the same.

Can mix use:

The legacy and new enclosure can mix use.

GE HealthCare

Technical Bulletin



New FRU:


Item	Legacy FRU No	New FRU No	Description	Qty	New FRU?
1	7500100	7500100-100	Cover Front RH Mulan Dual Port Locations-Asahi	1	Yes
2	5447195	5447195-100	COVER FRONT LH ASM 16 BEAT-Asahi	1	Yes
3	5444782-2	5444782-102	COVER FRONT TOP ASM 16BEAT-Asahi	1	Yes
4	5450848	5450848-100	Cover Rear LH Assembly 16 BEAT-Asahi	1	Yes
5	5441557	5441557-100	Cover Rear RH Assembly 16 BEAT-Asahi	1	Yes
6	5404661-2	5404661-102	COVER TOP, REAR, LH, BONDED, 16 BEAT-Asahi	1	Yes
7	5404663-2	5404663-202	COVER TOP, REAR , RH, BONDED, 16 BEAT-Asahi	1	Yes
8	5150411-2	5150411-102	Rear Trim Ring Assy for 16 BEAT-Asahi	1	Yes
9	5161159	5161159-100	Front Bridge Cover Assembly-Asahi	1	Yes
10	5404656-2	5404656-110	COVER SIDE LH BONDED,Signa High-Asahi. For 16 beat	2	Yes
11	5404657-2	5404657-110	COVER SIDE RH BONDED,Signa High-Asahi. For 16 beat	2	Yes
12	5457457	5457457-110	COVER SIDE CENTER RH ASM,Signa High-Asahi. For 16 Beat	1	Yes
13	5150410-2	5150410-102	Front Trim Ring Assy for 16 BEAT-Asahi	1	Yes
14	5404655	5404655-110	Side Cover, LH CTR, Bonded, Signa High-Asahi. For 16 Beat	1	Yes
15	5404658	5404658-100	COVER, TOP FRONT ASM 450w-Asahi	1	Yes

6. Affected Service Publications

A. SIGNA MR355/MR360 Service Methods (M:5856347-3EN R: 7/C:5856355-3EN R: 8/A:5856356-3EN R: 7)

Location: Home > Parts > Illustrated Parts > Enclosure and Gradient Assembly

Table 1. Enclosure and Gradient Assembly change to: **Yellow cells are new added.**

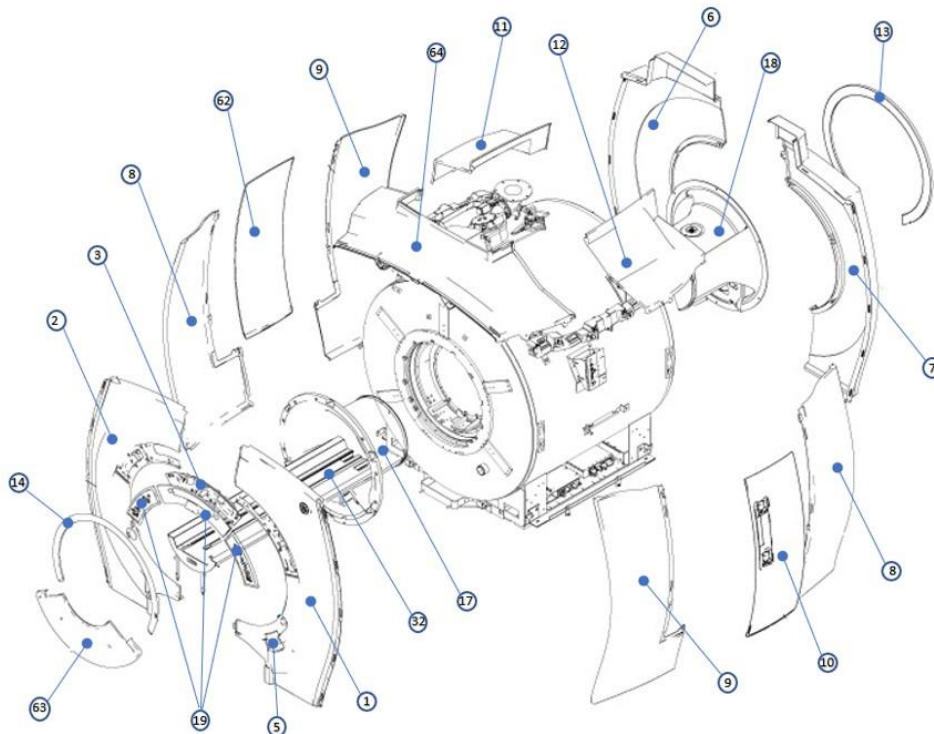

GE HealthCare

Technical Bulletin



Item	Part Number	Description	Qty	FRU	Lower
1	5446398	Front Cover Assmby RH	1	Yes	
	7500100	Cover Front RH Mulan Dual Port Locations	1	Yes	
	7500100-100	Cover Front RH Mulan Dual Port Locations-Asahi	1	Yes	
2	5447195	Front Cover Assmby LH	1	Yes	
	5447195-100	COVER FRONT LH ASM 16 BEAT-Asahi	1	Yes	
3	5444782-2	COVER FRONT TOP ASM 16BEAT	1	Yes	
	5444782-102	COVER FRONT TOP ASM 16BEAT-Asahi	1	Yes	
6	5450848	Rear cover assembly LH	1	Yes	
	5450848-100	Cover Rear LH Assembly 16 BEAT-Asahi	1	Yes	
7	5441557	Rear cover assm RH	1	Yes	
	5441557-100	Cover Rear RH Assembly 16 BEAT-Asahi	1	Yes	
8	5404656-2	Cover, Side Front bonded LH	2	Yes	
	5404656-110	COVER SIDE LH BONDED,Signa High-Asahi. For 16 beat	2	Yes	
9	5404657-2	Cover, Side Front bonded RH	2	Yes	
	5404657-110	COVER SIDE RH BONDED,Signa High-Asahi. For 16 beat	2	Yes	
10	5457457	Cover, Side center RH	1	Yes	
	5457457-110	COVER SIDE CENTER RH ASM,Signa High-Asahi. For 16 Beat	1	Yes	
11	5404661-2	Cover, Top rear LH ASSM	1	Yes	
	5404661-102	COVER TOP, REAR, LH, BONDED, 16 BEAT-Asahi	1	Yes	
12	5404663-2	Cover, Top rear RH ASSM	1	Yes	
	5404663-202	COVER TOP, REAR, RH, BONDED, 16 BEAT-Asahi	1	Yes	
13	5150411-2	Rear trim ring assm	1	Yes	
	5150411-102	Rear Trim Ring Assy for 16 BEAT-Asahi	1	Yes	
14	5150410-2	Front trim ring assm	1	Yes	
	5150410-102	Front Trim Ring Assy for 16 BEAT-Asahi	1	Yes	
15	5458261	Enclosure Mounting kit	1	Yes	
16	5180636	Turtle Bracket Asm Front BRM 1.5T	1	No	Turtle Bracket Asm Front BRM 1.5T
17	5160993-100	RoHS, Front End Bell BRM Assembly	1	Yes	Front End Bell BRM Assembly
18	5436428-100	RoHS, Rear End Bell Assy for 16 BEAT	1	Yes	Rear End Bell BRM Assembly
19	5447990	16 BEAT Control and Display Panel	1	No	Display / Control Panel
20	5142289-2	Gradient Cable Lead Board Molded Ring	1	Yes	
21	5142288-2	Gradient Cable Lead Board Cover	1	Yes	
22	2181231	FILLER PAD CONQUEST	4	Yes	
23	2344509-2	Bore Light End Cap, N9	2	Yes	
24	2222549	WIDE OPEN LUMITEX LIGHT PANEL 726 MM IN LENGTH	2	Yes	
25	2221945-100	CABLE PAC DISPLAY - ROHS	1	Yes	
26	2221948-51	Cable - SRI and display - RoHS	1	Yes	
27	5123813	Run#1261 - Operator Display Cable to SRI	1	Yes	
28	2221946	(862/863) MG2-A4-J1/J2 TO MG2-A5-J1/ J2	2	Yes	
29	2293539	TWIN BRIDGE SUPPORT BUSHING	4	Yes	
30	2293538-5	Split Bridge Edge Foam, Front, N9	2	Yes	
31	2226103	HOSE	2	Yes	

GE HealthCare

Technical Bulletin



32	5333181	SV front bridge assy	1	Yes	SV front bridge assy
33	2293202	TWIN BRIDGE SUPPORT FRT	1	Yes	
34	2293540	TWIN SUPPORT LOCK TOP	2	Yes	
35	2293541	TWIN BRIDGE SUPPORT SLIDER	2	Yes	
36	2259837	LCC Coldhead Control Cable Assembly	1	Yes	CABLE ASM,COLD HEAD CONTROL,LCC
37	2143771-2	2~ TEFLON TAPE	3	Yes	
38	2293771	RADIAL SHIM SPACER 1.5MM THK.	4	Yes	
39	2293771-2	RADIAL SHIM SPACER 0.8MM THK.	4	Yes	
40	5160992	Rear Bracket Assembly, Left	1	Yes	
41	5160992-2	Rear Bracket Assembly, Right	1	Yes	
42	46-208765P28	DIA RANGE 1-13/16 - 2-3/4, 302/305 SST ALL PARTS BREEZE CORP 30036 OR EQ	1	Yes	
43	46-170212P2	ADHESIVE THERMAL JOINT COMPOUND, 8-OZ JAR OPAQUE WHITE GREASE	0.01	Yes	
44	2185175	FOAM, CLOSED CELL, BRM CONQUEST	4	Yes	
45	46-294167P28	TUBE MENDER,CAT #91355K85, BRASS	2	Yes	
46	2174929-6	OPEN CELL POLYURETHANE FOAM, 45~ LG	3	Yes	
47	5118475-2	1.5T BRM-F Gradient Coil, GEHC N9	1	No	1.5T BRM-F Gradient Coil, GEHC N9
48	5432162	16 BEAT Platform Asm,PAC5 and SRI5 Mounting	1	No	Platform Asm, PAC and SRI Mounting
49	2379255	FOAM, GRAY TEXTURED SURFACE, 116~LG.	1	Yes	
50	2213168-3	Front Lower Support, BRM	1	Yes	
51	2213167-3	Rear Lower Support, BRM	1	Yes	
52	2240708	SUPPORT MOUNTING	2	Yes	
53	2238907	SPACER	1	Yes	
54	2238908	SPACER SHIM 1.5MM THK	2	Yes	
55	2238908-2	SPACER SHIM 2.4MM THK	6	Yes	
56	2241176	RUBBER ISOLATION PAD	2	Yes	
57	2238907-2	SPACER	1	Yes	
58	2212715-3	LCC BRM Mounting Pad	4	Yes	
59	5428034	Magnet Side Kit SV16CH	1	No	Magnet side KIT
60	5224349	Enclosure Cable Kit	1	No	Enclosure Cable Kit
61	2128699	(606) MAGNET END TO ERU/MRU		Yes	
62	5404655	Side Center LH	1	Yes	
	5404655-110	Side Cover, LH CTR, Bonded, Signa High-Asahi. For 16 Beat	1	Yes	
63	5161159	Bridge Cover	1	Yes	
	5161159-100	Front Bridge Cover Assembly-Asahi	1	Yes	
64	5404658	TOP FRONT	1	Yes	
	5404658-100	COVER, TOP FRONT ASM 450w-Asahi	1	Yes	
65	6134901	Name plate for MR380	1	Yes	

GE HealthCare


Technical Bulletin



B. Optima MR360/Brivo MR355/Optima MR360 Advance/Brivo MR355 Inspire Service Methods(SV23.0 platform & SV23.1 platform) (M:6336165-3EN R: 4, C:6336163-3EN R: 4, A:6336161-3EN R: 4)

Location: Content > Parts > Illustrated Parts > Enclosure and Gradient Assembly

Table 1. Enclosure and Gradient Assembly change to: **Yellow cells are new added.**

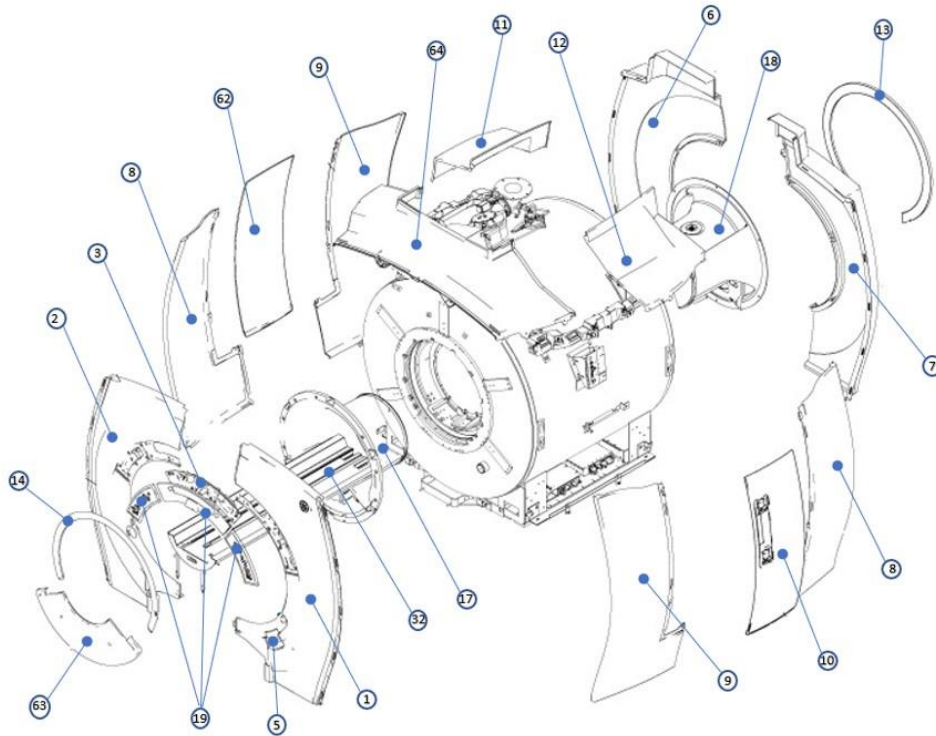

Item	Part Number	Description	Qty	FRU	Lower
1	5446398	Front Cover Assmbly RH	1	Yes	
	7500100	Cover Front RH Mulan Dual Port Locations	1	Yes	
	7500100-100	Cover Front RH Mulan Dual Port Locations-Asahi	1	Yes	
2	5447195	Front Cover Assmbly LH	1	Yes	
	5447195-100	COVER FRONT LH ASM 16 BEAT-Asahi	1	Yes	
3	5444782-2	COVER FRONT TOP ASM 16BEAT	1	Yes	
	5444782-102	COVER FRONT TOP ASM 16BEAT-Asahi	1	Yes	
6	5450848	Rear cover assembly LH	1	Yes	
	5450848-100	Cover Rear LH Assembly 16 BEAT-Asahi	1	Yes	
7	5441557	Rear cover assm RH	1	Yes	
	5441557-100	Cover Rear RH Assembly 16 BEAT-Asahi	1	Yes	
8	5404656-2	Cover, Side Front bonded LH	2	Yes	
	5404656-110	COVER SIDE LH BONDED,Signa High-Asahi. For 16 beat	2	Yes	
9	5404657-2	Cover, Side Front bonded RH	2	Yes	
	5404657-110	COVER SIDE RH BONDED,Signa High-Asahi. For 16 beat	2	Yes	
10	5457457	Cover, Side center RH	1	Yes	
	5457457-110	COVER SIDE CENTER RH ASM,Signa High-Asahi. For 16 Beat	1	Yes	
11	5404661-2	Cover, Top rear LH ASSM	1	Yes	

GE HealthCare

Technical Bulletin



	5404661-102	COVER TOP, REAR, LH, BONDED, 16 BEAT-Asahi	1	Yes	
12	5404663-2	Cover, Top rear RH ASSM	1	Yes	
	5404663-202	COVER TOP, REAR, RH, BONDED, 16 BEAT-Asahi	1	Yes	
13	5150411-2	Rear trim ring assm	1	Yes	
	5150411-102	Rear Trim Ring Assy for 16 BEAT-Asahi	1	Yes	
14	5150410-2	Front trim ring assm	1	Yes	
	5150410-102	Front Trim Ring Assy for 16 BEAT-Asahi	1	Yes	
15	5458261	Enclosure Mounting kit	1	Yes	
16	5180636	Turtle Bracket Asm Front BRM 1.5T	1	No	Turtle Bracket Asm Front BRM 1.5T
17	5160993-100	RoHS, Front End Bell BRM Assembly	1	Yes	Front End Bell BRM Assembly
18	5436428-100	RoHS, Rear End Bell Assy for 16 BEAT	1	Yes	Rear End Bell BRM Assembly
19	5447990	16 BEAT Control and Display Panel	1	No	Display / Control Panel
20	5142289-2	Gradient Cable Lead Board Molded Ring	1	Yes	
21	5142288-2	Gradient Cable Lead Board Cover	1	Yes	
22	2181231	FILLER PAD CONQUEST	4	Yes	
23	2344509-2	Bore Light End Cap, N9	2	Yes	
24	2222549	WIDE OPEN LUMITEX LIGHT PANEL 726 MM IN LENGTH	2	Yes	
25	2221945-100	CABLE PAC DISPLAY - ROHS	1	Yes	
26	2221948-51	Cable - SRI and display - RoHS	1	Yes	
27	5123813	Run#1261 - Operator Display Cable to SRI	1	Yes	
28	2221946	(862/863) MG2-A4-J1/J2 TO MG2-A5-J1/ J2	2	Yes	
29	2293539	TWIN BRIDGE SUPPORT BUSHING	4	Yes	
30	2293538-5	Split Bridge Edge Foam, Front, N9	2	Yes	
31	2226103	HOSE	2	Yes	
32	5333181	SV front bridge assy	1	Yes	SV front bridge assy
33	2293202	TWIN BRIDGE SUPPORT FRT	1	Yes	
34	2293540	TWIN SUPPORT LOCK TOP	2	Yes	
35	2293541	TWIN BRIDGE SUPPORT SLIDER	2	Yes	
36	2259837	LCC Coldhead Control Cable Assembly	1	Yes	CABLE ASM.COLD HEAD CONTROL LCC
37	2143771-2	2~ TEFLON TAPE	3	Yes	
38	2293771	RADIAL SHIM SPACER 1.5MM THK.	4	Yes	
39	2293771-2	RADIAL SHIM SPACER 0.8MM THK.	4	Yes	
40	5160992	Rear Bracket Assembly, Left	1	Yes	
41	5160992-2	Rear Bracket Assembly, Right	1	Yes	
42	46-208765P28	DIA RANGE 1-13/16 - 2-3/4, 302/305 SST ALL PARTS BREEZE CORP 30036 OR EQ	1	Yes	
43	46-170212P2	ADHESIVE THERMAL JOINT COMPOUND, 8-OZ JAR OPAQUE WHITE GREASE	0.01	Yes	
44	2185175	FOAM, CLOSED CELL, BRM CONQUEST	4	Yes	
45	46-294167P28	TUBE MENDER,CAT #91355K85, BRASS	2	Yes	
46	2174929-6	OPEN CELL POLYURETHANE FOAM, 45~ LG	3	Yes	
47	5118475-2	1.5T BRM-F Gradient Coil, GEHC N9	1	No	1.5T BRM-F Gradient Coil, GEHC N9
48	5432162	16 BEAT Platform Asm,PAC5 and SRI5 Mounting	1	No	Platform Asm, PAC and SRI Mounting
49	2379255	FOAM, GRAY TEXTURED SURFACE, 116~LG.	1	Yes	
50	2213168-3	Front Lower Support, BRM	1	Yes	

GE HealthCare

Technical Bulletin



51	2213167-3	Rear Lower Support, BRM	1	Yes	
52	2240708	SUPPORT MOUNTING	2	Yes	
53	2238907	SPACER	1	Yes	
54	2238908	SPACER SHIM 1.5MM THK	2	Yes	
55	2238908-2	SPACER SHIM 2.4MM THK	6	Yes	
56	2241176	RUBBER ISOLATION PAD	2	Yes	
57	2238907-2	SPACER	1	Yes	
58	2212715-3	LCC BRM Mounting Pad	4	Yes	
59	5428034	Magnet Side Kit SV16CH	1	No	Magnet side KIT
60	5224349	Enclosure Cable Kit	1	No	Enclosure Cable Kit
61	2128699	(606) MAGNET END TO ERU/MRU		Yes	
62	5404655	Side Center LH	1	Yes	
	5404655-110	Side Cover, LH CTR, Bonded, Signa High-Asahi. For 16 Beat	1	Yes	
63	5161159	Bridge Cover	1	Yes	
	5161159-100	Front Bridge Cover Assembly-Asahi	1	Yes	
64	5404658	TOP FRONT	1	Yes	
	5404658-100	COVER, TOP FRONT ASM 450w-Asahi	1	Yes	
65	6134901	Name plate for MR380	1	Yes	

Revision	Version	Date	Author	Description
01	01	15-March-2024	SONG Yang (SSO: 503386740)	Initial Release

End of Document