

MAGNETOM

Software System

Software Installation syngo MR B15 (NUMARIS/4 VB15A)

Document Version

Siemens reserves the right to change its products and services at any time.

In addition, manuals are subject to change without notice. The hardcopy documents correspond to the version at the time of system delivery and/or printout. Versions to hardcopy documentation are not automatically distributed.

Please contact your local Siemens office to order a current version or refer to our website <http://www.healthcare.siemens.com>.

Disclaimer

Siemens provides this documentation "as is" without the assumption of any liability under any theory of law.

The content described herein requires superior understanding of our equipment and may only be performed by qualified personnel who are specially trained for such installation and/or service.

Copyright

"© Siemens, 2007" refers to the copyright of a Siemens entity such as:

Siemens Healthcare GmbH - Germany

Siemens Aktiengesellschaft - Germany

Siemens Shenzhen Magnetic Resonance Ltd. - China

Siemens Shanghai Medical Equipment Ltd. - China

Siemens Healthcare Private Ltd. - India

Siemens Medical Solutions USA Inc. - USA

Siemens Healthcare Diagnostics Inc. - USA and/or

Siemens Healthcare Diagnostics Products GmbH - Germany

1	General	5
1.1	syngo MR system software	5
1.1.1	Introduction	5
1.1.2	Prerequisites	5
1.1.3	Abbreviations	5
2	Installing MRC / MRC+MRSC	6
2.1	Overview of installing software MRC/MRC+MRSC R630/640	6
2.2	Software kits	7
2.2.1	Celsius R630 4GB	7
2.2.2	Celsius R640 4GB	7
2.3	BIOS Update	8
2.3.1	Load BIOS defaults for Host	8
2.4	NUMARIS/4 Installation (host without operating system)	9
2.5	NUMARIS/4 reinstallation (on booting HW)	10
2.6	NUMARIS/4 software installation on R650	12
2.6.1	Overview of installing software MRC/MRC+MRSC	12
2.6.2	BIOS update	13
2.6.3	NUMARIS/4 installation (host without operating system)	13
2.6.4	NUMARIS/4 reinstallation (on booting HW)	14
2.7	Imager	16
2.7.1	BIOS Update	16
3	Site-specific data	17
3.1	Backup site-specific data on CD	17
3.2	Backing up site-specific data on external USB device	18
3.3	Restore site-specific data (manually after SW installation)	19
3.4	Restore site-specific from external USB device	20
3.4.1	Final checks	20
4	Configuration	21
4.1	General	21
5	Appendix	22
5.1	Language settings	22
5.1.1	Keyboard and Language setting	22
5.1.2	Language Settings Japanese	26
5.1.3	Language Settings Chinese	34
5.2	Archiving	44
5.2.1	Configuring the CD-R drive as an Archive Device	44
5.2.2	Settings for Codonix Color printer NP1660M	45
5.3	Configuration	46
5.3.1	Chorus	46
5.3.2	Chorus and PhoenixZIP	47

5.3.3	Enable/Disable Inline Image Filter	49
5.4	Troubleshooting	51
5.4.1	Customer network IP conflict with measurement network	51
5.4.2	Problems during syngo Installation	51
5.4.3	Workaround for Satellite Console connection problem	52
5.4.4	Problems with MRSC Start up	53
6	Changes	55
7	List of Hazard IDs	56

1.1 syngo MR system software

1.1.1 Introduction

The installation and service of equipment described herein is to be performed by qualified personnel employed or authorized by Siemens Medical Solutions or one of its affiliates. Assemblers and other personnel not employed by nor directly affiliated with nor authorized by Siemens Medical Solution or one of its affiliates` technical services are directed to contact the local office of Siemens Medical Solutions or its affiliates before attempting installation or service procedures.

1.1.2 Prerequisites

To configure the syngo-based MR system, the following training is required in combination with this software manual:

- syngo training / syngo MR training
- Network training (TCP/IP)
- DICOM training



All passwords required to install and configure the system are available in the "CS Knowledge Base" <https://smskb.medical.siemens.com/selfservice> under "Password List" (except for the service key (variable password), which is handled by the districts).

1.1.3 Abbreviations

EDB	Event Data Base
n.a.	Not applicable
SRS	Siemens Remote Service
QDB	Quality Assurance Data Base

2.1 Overview of installing software MRC/MRC+MRSC R630/640

Following table provides some important milestones for the software installation.



The software media is host type dependent, make sure the correct kit is used for your Host type (Celsius R630 or 640).

MRC (save set available)	MRC (no valid save set)	MRSC
Install software	Install software	Install software
Restore site specific data	Initial configuration	Restore site specific data / or import clone package from the MRC
Check System language and weight	Configure the system	Check System language and weight
Check if Service packs are available for this software version	Perform system tune-up, Quality Assurance, Installation Protocol	Check if Service packs are available for this software version
Check the image orientation	Check System language and weight	n.a.
Check if Update Instructions with service packs are available for this software version		
Perform the BIOS Update for the Imager		

NOTICE

Verify the software baseline label with the function Help, Info.

Correct software label for the released syngo MR B15 is:

- » COEM VE10B
- » syngo VE31G
- » syngo VX69A
- » N4_VB15A_LATEST_20070519
- » N4_VB15A_LATEST_20070519_P2
- » N4_VB15A_LATEST_20070519_P12
- » VX69A SL02P02

2.2 Software kits

2.2.1 Celsius R630 4GB

DVD-ROM (3 DVDs)

syngo MR B15

NUMARIS/4 VB15A

MAGNETOM Avanto / Espree / Symphony a Tim System / Trio a Tim System

Revision A

with part number **10277193** K2219

DVD1/3:Software R630

DVD2/3:Source DVD

DVD3/3:Service-pack 1

2.2.2 Celsius R640 4GB

DVD-ROM (3 DVDs)

syngo MR B15

NUMARIS/4 VB15A

MAGNETOM Avanto / Espree / Symphony a Tim System / Trio a Tim System

Revision A

with part number **10276988** K2219

DVD1/3:Software R640

DVD2/3:Source DVD

DVD3/3:Service-pack 1

2.3 BIOS Update

The BIOS update optimizes the setting for NUMARIS/4 for your host type, sets the NUMA-RIS/4 BIOS password and boot device order

2.3.1 Load BIOS defaults for Host

- Reboot the system (power on the Host).
- During the boot press **F2** to access the BIOS setup
- Enter the BIOS-Password.
- Check date / time, correct if necessary. Installations in Asia should change the date 3 days back before the installation in order to avoid problems with the licenses.
- Check in Security that entry Start from IDE CD-ROM is set to Enabled.
- Select the **Boot**-Tab.
- Highlight the **Boot Device Priority** and press **Enter**.
- Move the CD-ROM-Entry to the top of the list.
- If there is an ! next to the entry, then press SHIFT+ Insert BIOS-CD VB15x: xxxxxxxx-MRxxxBIOS-CD in to the DVD/CD-ROM drive.
- Press **F10** and accept "Yes" by pressing **Enter**.
- The BIOS-update will run, wait until the message "Ready for Software Installation" is displayed.
- Remove the BIOS-CD after the automatic eject.
- Restart the system by pressing **Ctrl + Alt + Del** simultaneously.

2.4 NUMARIS/4 Installation (host without operating system)

Following table describes how to install software on hardware without operating system.

- Perform the "BIOS update"

Action	MRC	MRSC
Insert the "Numaris/4 DVD 1/3" in to the DVD/CD-ROM (R:) drive.		
Insert the "Master-Backup-CD" to the CD/DVD-burner drive DVD/CD-burner (S:)		
... automatic boot will follow		
A pseudo-graphic user interface is displayed with the warning that all data will be lost after pressing a key for confirmation. At this point, press Enter to accept the warning that all data will be destroyed.		
... automatic reboot will follow		
A pseudo-graphic menu with the following windows will be displayed: "Selection Menu", "Hardware Detection Results", and "Information Window"		
Press Enter twice (2 times) to continue the installation.		
... several reboots will follow		
... Symantec Ghost window will be displayed with a "Progress Indicator" this will take ~25 min.		
Remove the "syngoMR xxx Customer Image" DVD from the CD/DVD-ROM drive (R:). When the following message is displayed: Plug and Play is running, please wait! Depending on the hardware, this may take 2-5 minutes		
Insert the Service Pack CD into the CD/DVD-ROM drive (R:)		
The automatic installation will continue, this will take ~25 min. At the end of the installation the "Windows XP Login window" will be displayed.		
Log in as "meduser"		
Click on the link Click here to start the Local Service Configuration in the Internet Explorer (with the window title Start Up Check).		
Enter the valid Service Key and Password, click OK .		
Select Configuration , navigate through the configuration, check that the network settings and host name are correct.		

2.5 NUMARIS/4 reinstallation (on booting HW)

Following table describes how to install software on hardware that is bootable.

- Switch on the Host

Action	MRC	MRSC
Login as administrator		
Insert the "Numaris/4 ... DVD 1/3" in to the DVD/CD-ROM drive.		
Insert the "Master-Backup-CD" to the CD/DVD-burner drive (DVD/CD-burner)		
Navigate to "R:\start.bat" (use the file explorer)		
Double-click "start.bat"		
A Dos-shell pops up with the following message: Start to irreversible prepare the system for installation [Y/N]?		
Press Y to continue		
... automatic reboot will follow		
A pseudo-graphic user interface is displayed with the warning that all data will be lost after pressing a key for confirmation. At this point, press Enter to accept the warning that all data will be destroyed.		
... automatic reboot will follow		
A pseudo-graphic menu with the following windows will be displayed: "Selection Menu", "Hardware Detection Results", and "Information Window"		
Press Enter twice (2 times) to continue the installation.		
... several reboots will follow		
... Symantec Ghost window will be displayed with a "Progress Indicator" this will take ~25 min.		
Remove the "syngoMR xxx Customer Image" DVD from the CD/DVD-ROM drive (R:\). When the following message is displayed: Plug and Play is running, please wait! Depending on the hardware, this may take 2-5 minutes		
Insert the Service Pack CD into the CD/DVD-ROM drive (R:\)		
The automatic installation will continue, this will take ~25 min. At the end of the installation the "Windows XP Login window" will be displayed.		
Log in as "meduser"		
Click on the link Click here to start the Local Service Configuration in the Internet Explorer (with the window title Start Up Check).		

Action	MRC	MRSC
Enter the valid Service Key and Password, click OK .		
Select Configuration , navigate through the configuration, check that the network settings and host name are correct.		

2.6 NUMARIS/4 software installation on R650

2.6.1 Overview of installing software MRC/MRC+MRSC

The following table provides a number of important milestones for the software installation.

MRC (save set available)	MRC (no valid save set)	MRSC
Install software	Install software	Install software
Install software at the Imager		n.a.
Restore site-specific data	Initial configuration	Restore site-specific data / or import clone package from the MRC
Check system language and weight	Configure the system	Check system language and weight
Check if service packs are available for this software version	Perform system tune-up, quality assurance, and the installation protocol	Check if service packs are available for this software version
Check the image orientation	Check system language and weight	n.a.
Check if Update Instructions with service packs are available for this software version		

NOTICE

Verify the software baseline label with the function Help, Info.

Correct software label for the released syngo MR B15 is:

- » COEM VE10B
- » syngo VE31G
- » syngo VX69A
- » Revision B
- » N4_VB15A_LATEST_20070519
- » N4_VB15A_LATEST_20070519_P2
- » N4_VB15A_LATEST_20070519_P12
- » VX69A SL02P02

2.6.2 BIOS update

The BIOS update optimizes the setting of NUMARIS/4 for your host type, sets the NUMARIS/4 BIOS password and boot device order

2.6.2.1 Loading BIOS defaults for Celsius R650

- Reboot the system (power on the Host).
- During the boot, press **F2** to access the BIOS setup
- Enter the BIOS password.
- Check date / time, correct if necessary. Installations in Asia should change the date going back 3 days for the installation to avoid problems with the licenses.
- Check in Security that the entry Start from CD-ROM is set to Enabled.
- Select the **Boot** tab.
- Highlight the **Boot Device Priority** and press **Enter**.
- Move the CD-ROM entry to the top of the list.
- Press **F9** and accept (Yes) by pressing **Enter**.
- If you see an ! next to the entry, press SHIFT+ Insert BIOS-CD Vxxx: xxxxxxx-MR005BIOS-CD in to the DVD/CD-ROM drive.
- Press **F10** and accept "Yes" by pressing **Enter**.
- The BIOS update will run, wait until the "Ready for Software Installation" message is displayed.
- Remove the BIOS CD after the automatic eject.
- Restart the system by pressing **Ctrl + Alt + Del** simultaneously.

2.6.3 NUMARIS/4 installation (host without operating system)

The following table describes how to install software on hardware without an operating system.

- Perform the "BIOS update"

Action	MRC	MRSC
Insert the "Numaris/4 DVD 1/2" in to the DVD/CD-ROM (R:) drive.		
Restart the system by pressing Ctrl + Alt + Del simultaneously.		
... automatic boot will follow		
The installation environment WINPE is booted from the DVD.		
Wait until you get the Destroy data popup		
Confirm the displayed popup Destroy Data on the partition C: with Yes		
Now the partitions are automatically destroyed created new.		

Action	MRC	MRSC
... Symantec ghost window will be displayed with a "progress indicator" and all partition are restored from the DVD image. The overall procedure will take ~50 min. After the sucessfull restore of the images a reboot is performed automatical		
Remove the "syngoMR xxx Customer Image" DVD from the CD/DVD-ROM drive (R:\). When the following message is displayed: Plug and Play is running, please wait! Depending on the hardware, this may take 2-5 minutes		
Insert the "Master Backup-CD" into the CD/DVD burner drive DVD/CD burner (S:)		
Insert the Service Pack CD into the CD/DVD-ROM drive (R:\)		
The automatic installation will continue, this will take ~25 min. At the end of the installation, the "Windows XP login window" will be displayed.		
Log in as "meduser"		
Click the link Click here to start the Local Service Configuration in the Internet Explorer (with the window title Start Up Check).		
Enter the valid Service Key and Password, click OK .		
Select Configuration , navigate through the configuration, check to ensure that the network settings and host name are correct.		

2.6.4 NUMARIS/4 reinstallation (on booting HW)



For the software reinstallation the backup cd should be not in the CD/DVD-R drive [S] during the software installation. If the backup cd is in CD/DVD-R drive [S] the software installation from the DVD install media fails.

You have to follow the procedure shown below.



Ignore the message "Zlib Exception" during the restore of the backup data.

The following table describes how to install software on hardware that is bootable.

- Switch on the Host

Action	MRC	MRSC
Login as administrator		
Insert the "Numaris/4 ... DVD 1/2" into the DVD/CD-ROM drive.		

Action	MRC	MRSC
Navigate to "R:\start.bat" (use File Explorer)		
Double-click "start.bat"		
... automatic boot will follow		
The installation environment WINPE is booted from the DVD.		
Wait until you get the Destroy data popup		
Confirm the displayed popup Destroy Data on the partition C: with Yes		
Now the partitions are automatically destroyed created new.		
... Symantec ghost window will be displayed with a "progress indicator" and all partition are restored from the DVD image. The overall procedure will take ~50 min. After the successful restore of the images a reboot is performed automatic		
Remove the "syngoMR xxx Customer Image" DVD from the CD/DVD-ROM drive (R:). When the following message is displayed: Plug and Play is running, please wait! Depending on the hardware, this may take 2-5 minutes		
Insert the "Master Backup-CD" into the CD/DVD burner drive (DVD/CD burner)		
Insert the Service Pack CD into the CD/DVD-ROM drive (R:)		
The automatic installation will continue, this will take ~25 min. At the end of the installation, the "Windows XP login window" will be displayed.		
Log in as "meduser"		
Click the link Click here to start the Local Service Configuration in the Internet Explorer (with the window title Start Up Check).		
Enter the valid Service Key and Password, click OK .		
Select Configuration , navigate through the configuration, check to ensure that the network settings and host name are correct.		

2.7 Imager

2.7.1 BIOS Update



Before performing the BIOS update at the Imager, make sure that the most recent service pack is installed at the MRC.

Action	MRC
Login as administrator	
Make sure that the Imager is connected to the MRC	
Insert the BIOS-CD VB15A-BIOS-CD into DVD-ROM Drive.	
Navigate to "R:\BootNet\Host" (use the file explorer)	
Double-click "BootNet.bat"	
Restart the imager as demanded in the pop-up-message.	
Wait for the Actualisation-process to complete	
Reboot the imager as demanded.	
Remove the BIOS-CD from the MRC	



If the BIOS-CD VB15A is used for a second time to update the Imager, only then perform the following steps after the installation:

Action (only valid for 2:nd installation of Imager BIOS)	MRC
On the MRC navigate to C:\MRIR\Med\Config\scripts	
Delete the file File: setup.lastcall.	
Reboot the Imager	

3.1 Backup site-specific data on CD

Required time: 10-60 min

Required tools / parts: CD-R



Use only CD-R 's with medical grade, e.g. Siemens Medical CD-R 80 or TDK Medical CD-R.

Save the site-specific data each time you make changes to the site-specific configuration, tune-up, or sequences and protocols. In case of a hard disk crash, you will be able to restore a valid set of site-specific data.

- Login as **administrator**
- Start the service platform, select **Backup & Restore**.
- Insert a recordable CD-disk in the CD/DVD-burner.
- Select **Backup** in the "Command" pull-down menu.
- In the menu Drives, select the target drive CD/DVD-burner (drive letter **S:**) used to save the data.
- Select **Master-Backup** in the "Packages" menu.
- Click **Go** to start the "Backup", check for errors in the log-window after the backup has been completed.

Setting that are not part of the site-specific data:

Remote access mode (No Access, Limited Access, etc. in the customer UI).

Auto-transfer flag ("Configuration", "Service", "Auto Transfers")

Postscript Printer Configuration (Customer paper printer)

ImageText Settings (ImageText Settings from ImageText Dialog)

Time zone

Transfer Rules



During backup at the MRSC, a "**CorRea.bru**" error may be displayed, please ignore this error.

3.2 Backing up site-specific data on external USB device



Please check the external media (e.g. USB-disk) first with the virus scanner of the service PC. Be sure that the external media used is free of viruses.

- Login as **administrator** ↵.
- Start the "Internet Explorer"
- Connect the USB device.
- Check what drive letter was set for the USB device with the windows file explorer.
- Login to the service platform, select **Backup & Restore**.
- Select Backup in the Command menu, select drive letter of the USB device e.g. **F** in the Drives menu..
- Select **Master-Backup** in the Package menu, start the backup by clicking **OK**.
- (Estimated time ~10-30 min) check for errors using the log window.
- Click **Home** to leave Backup & Restore.

Setting that are not part of the site-specific data:

Remote access mode (No Access, Limited Access, etc. in the customer UI).

Auto-transfer flag ("Configuration", "Service", "Auto Transfers")

Postscript Printer Configuration (Customer paper printer)

ImageText Settings (ImageText Settings from ImageText Dialog)

Time zone

Transfer Rules



During backup at the MRSC, a "**CorRea.bru**" error may be displayed, please ignore this error.

3.3 Restore site-specific data (manually after SW installation)

Restoring the site-specific data may be helpful in many situations. Please read the following information before beginning the work.

After a software installation, a set of site-specific data will save you plenty of work/time. The site-specific data e.g. Configuration, Tune-up, and QA are included in the master backup.

- Login as **meduser** ↵.
- Start the "Internet Explorer"
- Click the link "Click here to start the Local Service Configuration" in the Internet Explorer (with the window title Start Up Check).
- Login to the service platform, select **Backup & Restore**.
- Insert the backup disk in the CD/DVD-ROM.
- Select **Restore** in the "Command" pull-down menu.
- Select the source drive CD/DVD-ROM (drive letter **R:**) in the "Drives" menu.
- Select your xxx_**Master-Backup**_xxxx_xxxx_ok.ar to be restored in the "Archive" menu.
- In the menu "Groups", select the necessary packages to be restored.
- Click **Go** to start the "Backup", check for errors in the log window after restore has been completed. (Estimated time ~10-15 min.)
- Click **Home**.
- Click **OK** to reboot the system.



For syngo MR Espree systems the "Tune-up" step "Coil Power Losses" must be performed after restore of site-specific data.



To ensure the Administrator account password is correct following restoration of the dynamic data, it is recommended to set the default Administrator password in the Service GUI > **Configuration** > **Local Host** > **Users**

3.4 Restore site-specific from external USB device



Please check the external media (e.g. USB-disk) first with the virus scanner of the service PC. Be sure that the external media used is free of viruses.

Restoring the site-specific data may be helpful in many situations. Please read the following information before beginning with the work.

In case of corrupted site-specific data, we recommend saving the actual status of the site-specific data to another medium before restoring an old save set. This will save time if the problem is resolved by reinstalling the software. The actual tune-up data may save a great deal of work, especially if the customer optimized his protocols recently.

- Login as **administrator** ↵.
- Start the "Internet Explorer"
- Login to the service platform, select **Backup & Restore**.
- Connect the USB device.
- Select **Restore** in the "Command" drop-down menu.
- Select the source USB-drive (drive letter e.g. **F:**) in the "Drives" menu.
- Select your xxx_**Master-Backup**_xxxx_xxxx_ok.ar to be restored in the "Archive" menu.
- In the menu "Groups", select the necessary packages to be restored.
- Click **Go** to start the "Backup", check for errors in the log window after restore has been completed. (Estimated time ~20-60 min.)



For syngo MR Espree systems the "Tune-up" step "Coil Power Losses" must be performed after restore of site-specific data.



To ensure the Administrator account password is correct following restoration of the dynamic data, it is recommended to set the default Administrator password in the Service GUI > **Configuration** > **Local Host** > **Users**

3.4.1 Final checks

- Log in as **meduser** ↵.
- Perform the "Quality Assurance" procedure: "Image Orientation" in the service platform..
- Make sure that the system is generating images.

4.1 General

For detailed information regarding the configuration of syngo MR systems, refer to **Help** in the Service Software Platform

If you have changed the configuration within the Service Software Platform and the message "reboot system pending" is displayed in the "status bar", reboot the system as follows: Click **Home**, confirm the reboot message window with **OK**.



Language settings and the configuration UI for option Chorus are described in chapter Appendix.



How to configure Siemens Remote Service is described separately in the document "Installation and Startup for SRS".



The use of non-released combinations of hardcopy cameras and BU products is prohibited and does not meet the requirements defined. For detailed information regarding the CAMERA configuration, refer to technical documentation "Camera Info ..." valid for your NUMARIS/4 software version.



The "CCC" Connectivity Competence Center provides a powerful configuration management tool see:

<https://cs.med.siemens.de/cms/en/home/>

For Service, DICOM, in the DICOM Homepage (CCC) navigate to "Databases", "Interoperability DB", **Download, Acquisition Tool**).



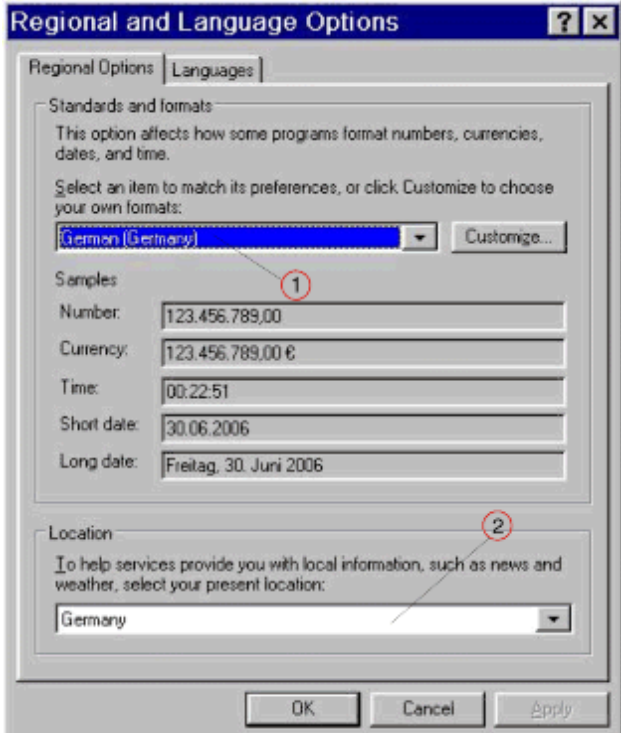
Remember to create a backup if the configuration has changed.

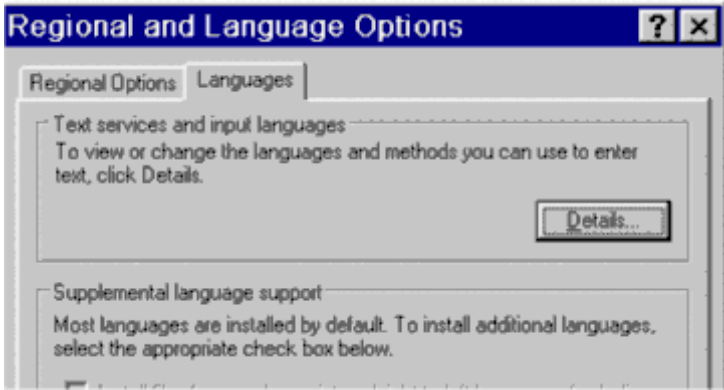
5.1 Language settings

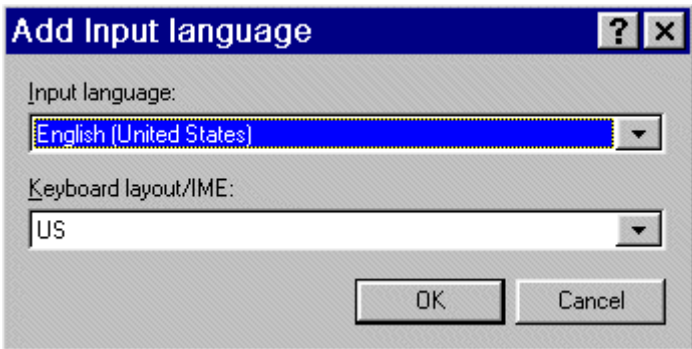
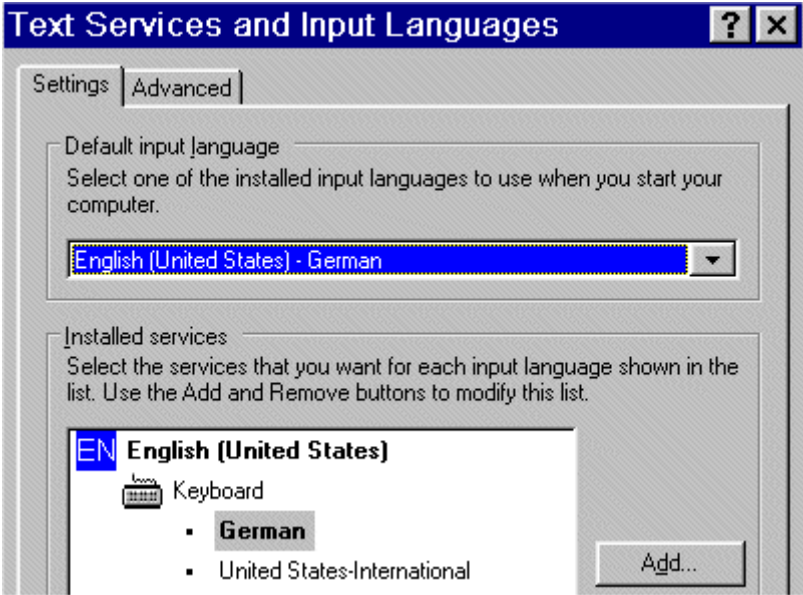
5.1.1 Keyboard and Language setting

- Setting the Keyboard Layout and Language

ACTION	MRC	MRSC
■ In the Main UI select Options -> Configuration	☐	☐

ACTION	MRC	MRSC																				
■ Click OK. Select the language for the main UI from the pull-down menu. The table below shows the supported entries for the languages. Under location, you can choose your country. Fig. 1:  (1) Language (2) Country Only the following entries can be used for 1 in the graphic above: <table><tr><th>main UI language</th><th>entry in the pull down menu</th></tr><tr><td>german</td><td>German (Germany)</td></tr><tr><td>english</td><td>English (US)</td></tr><tr><td>english (metric)</td><td>English (UK)</td></tr><tr><td>french</td><td>French (France)</td></tr><tr><td>spanish</td><td>Spanish (Spain)</td></tr><tr><td>chinese</td><td>Chinese (PRC)</td></tr><tr><td>japanese</td><td>Japanese</td></tr></table> <tr><td></td><td>☐</td><td>☐</td></tr> <tr><td>■ Click Apply.</td><td>☐</td><td>☐</td></tr>	main UI language	entry in the pull down menu	german	German (Germany)	english	English (US)	english (metric)	English (UK)	french	French (France)	spanish	Spanish (Spain)	chinese	Chinese (PRC)	japanese	Japanese		☐	☐	■ Click Apply .	☐	☐
main UI language	entry in the pull down menu																					
german	German (Germany)																					
english	English (US)																					
english (metric)	English (UK)																					
french	French (France)																					
spanish	Spanish (Spain)																					
chinese	Chinese (PRC)																					
japanese	Japanese																					
	☐	☐																				
■ Click Apply .	☐	☐																				

ACTION	MRC	MRSC
■ Select the Language tab in the Regional and Language Options task card. Fig. 2: 	☐	☐
■ Click Details.	☐	☐
■ In the upcoming window "Text Services and Input Languages" press Add Fig. 3: 	☐	☐

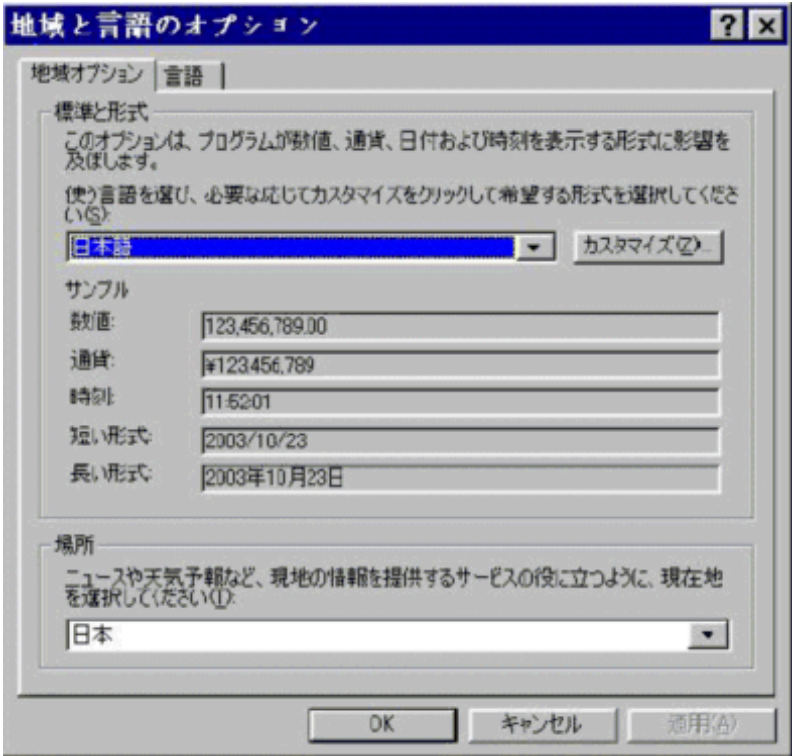
ACTION	MRC	MRSC
- In the next pop-up window choose the keyboard layout used under "Keyboard Layout." Select the language used for the Main UI under "Input language." The input language also determines the used system weight unit. Fig. 4:  - English (US) uses pounds as the unit of weight and inches as the unit of length. - English (UK) uses kg as the unit of weight and cm as the unit of height.	☐	☐
Click OK.	☐	☐
Click Apply.	☐	☐
In the "Text Services and Input Languages" window choose the newly configured keyboard layout under "default input language." Fig. 5: 	☐	☐
Click Apply	☐	☐

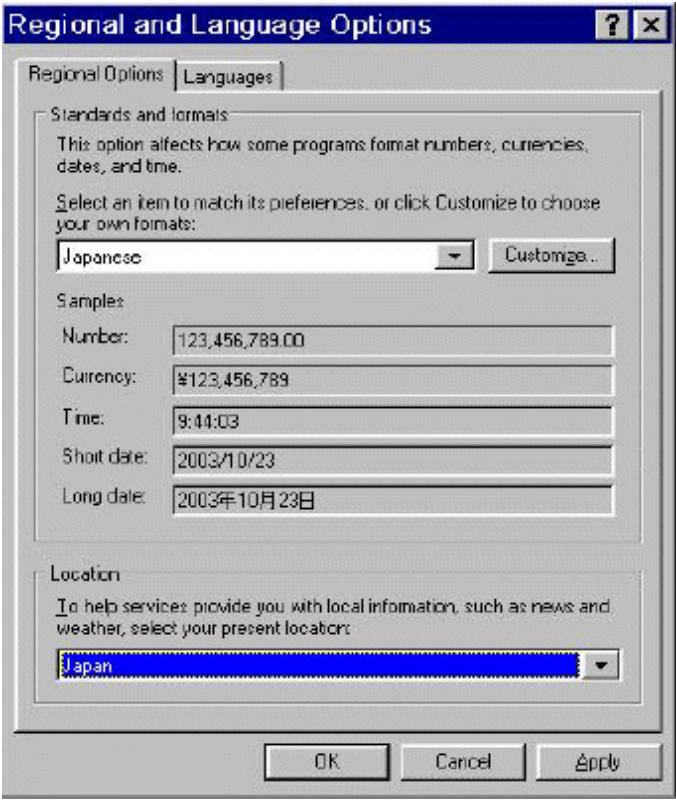
5.1.2 Language Settings Japanese

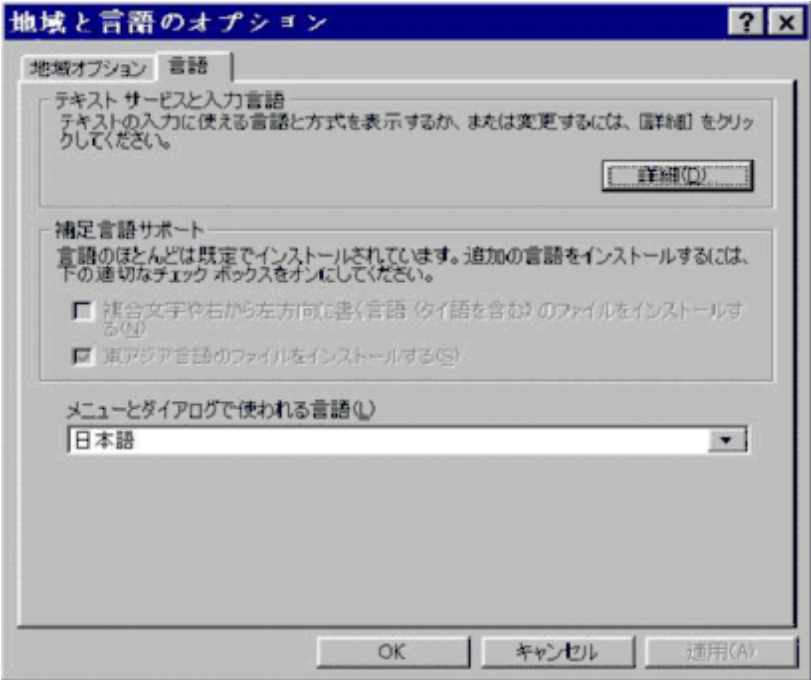
- Switch the language from international (US English) to Japanese

ACTION	MRC	MRSC
Select Options -> Configuration in the main menu of the Meduser platform	☐	☐
Double-click the Regional and Language Settings icon	☐	☐
In the highlighted field, choose Output pop-up in Japanese, if it was already switched during the Windows XP installation	☐	☐

Fig. 6:



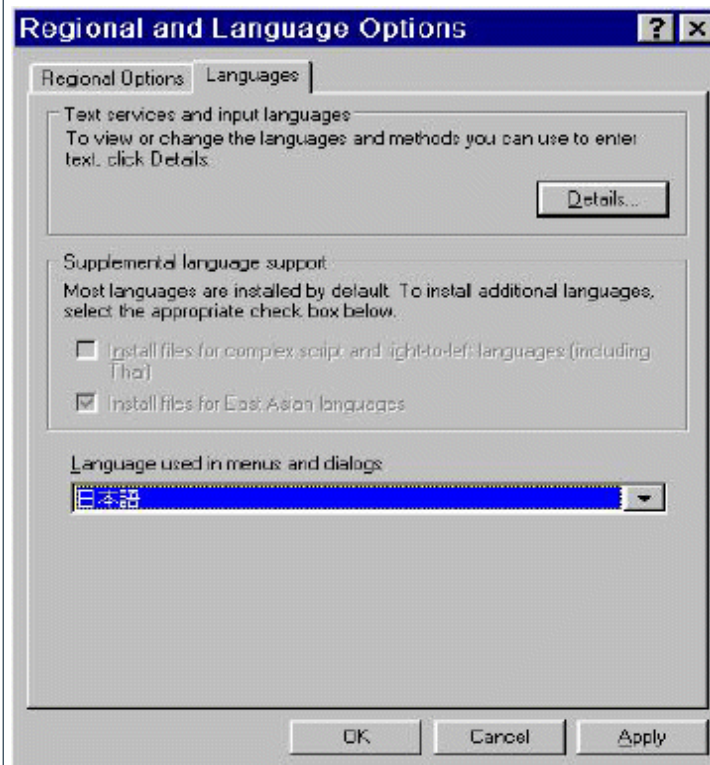
ACTION	MRC	MRSC
Otherwise, the international view of the pop-up is displayed. Fig. 7:  The first item (see below) Fig. 8:  (Remark: scroll up, if the following symbol is not present) Note: Be sure to use the long format for the date. It must look like yyyy/mm/dd or yyyy-mm-dd. If it is not set, you have to change it. For this purpose, press the Customize button and select the correct setting for the date.		

ACTION	MRC	MRSC
■ Switch to the 2. tab card "language" and choose the same Japanese symbol as in the first tab card Japanese symbol: Fig. 9:  Japanese view Fig. 10: 	☐	☐

ACTION

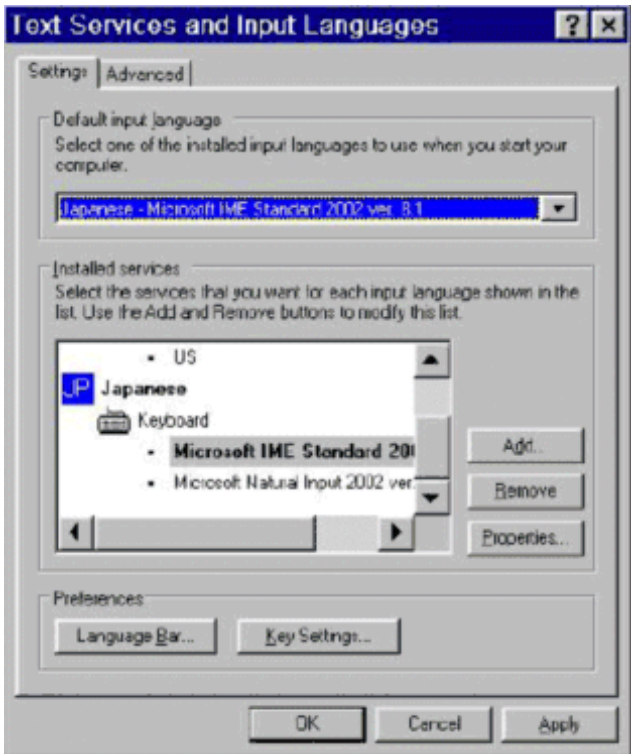
or international (US English) view.

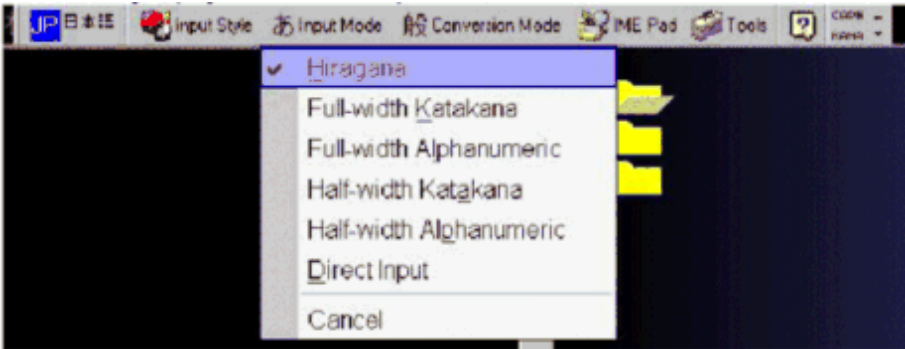
Fig. 11:



MRC

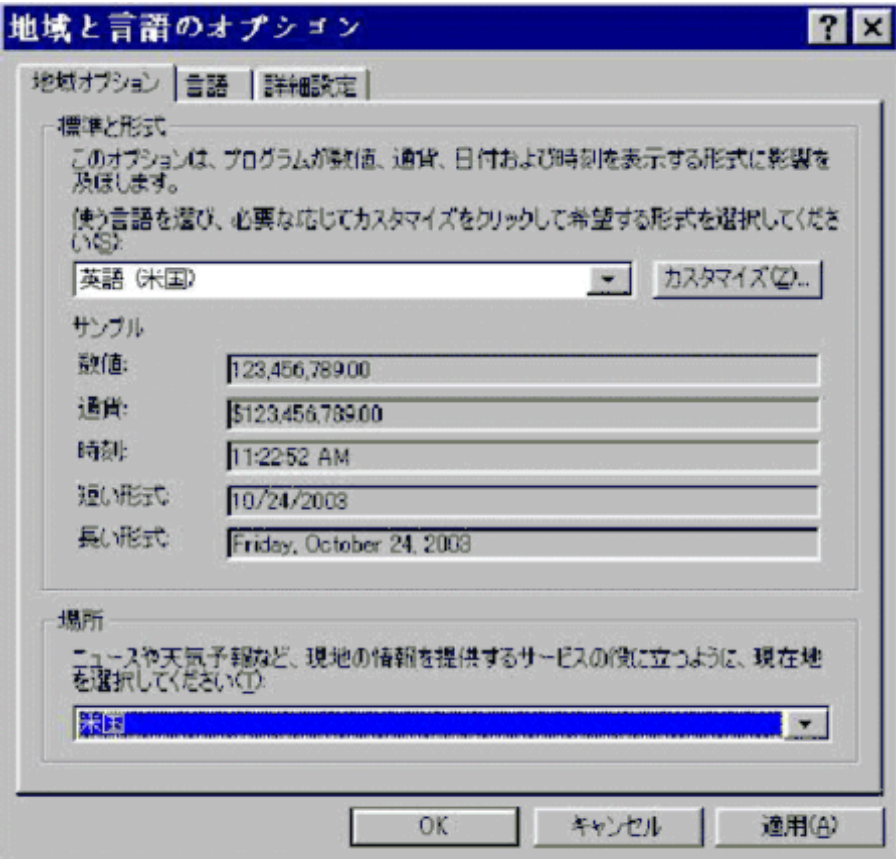
MRSC

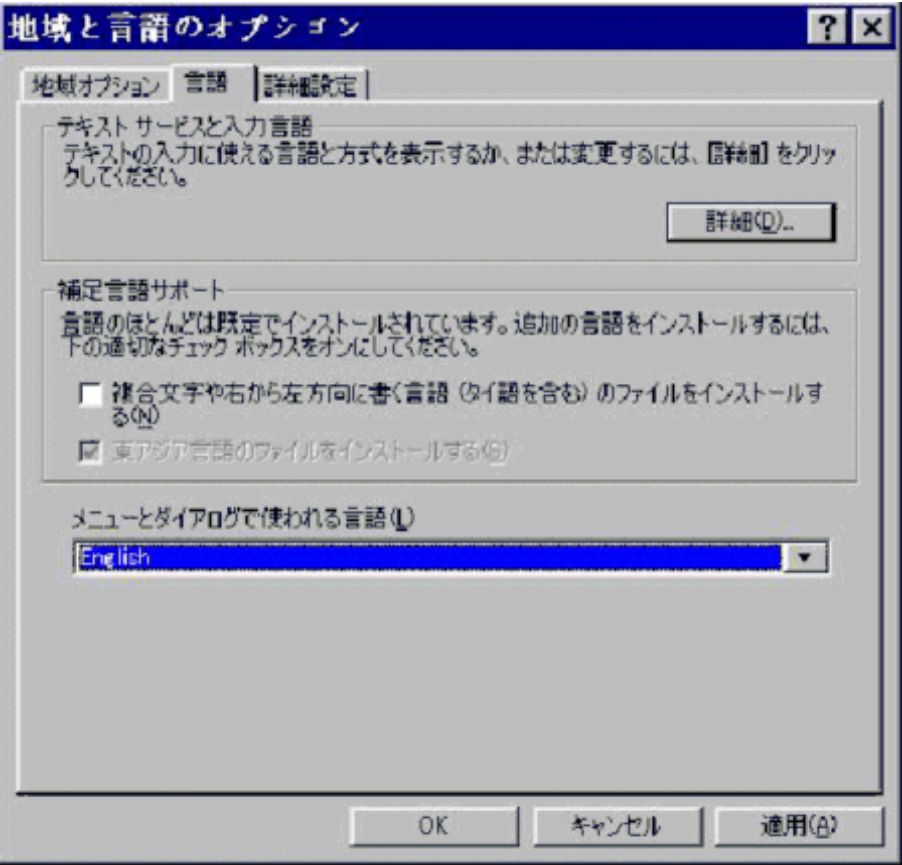
ACTION	MRC	MRSC
■ Click the "Details" button and perform the following settings for your keyboard: Japanese view Fig. 12:  or international (US English) view Fig. 13: 	☐	☐

ACTION	MRC	MRSC
IME settings: (right-click "Input Mode").	☐	☐
- Be sure that the Dicom Character Set is set to "ISO 2002 IR13" for Japanese in the Configuration (local service)! Fig. 16:  The screenshot shows a software interface with a menu titled 'Input Mode'. The menu is open, displaying several options: 'HIRAGANA' (which is highlighted with a blue bar and a checkmark), 'Full-width Katakana', 'Full-width Alphanumeric', 'Half-width Katakana', 'Half-width Alphanumeric', 'Direct Input', and 'Cancel'. The background of the software interface is dark, and there are various icons and text elements visible around the menu. - Start local service - Go in Configuration - Select the link 'Character Set' under DICOM - Choose the radio button "allow multiple character sets" - Select the Character Set "ISO 2002 IR13" - Press "Save" and confirm the pop-up. - Restart the system by pressing "Home." Confirm the restart pop-ups.	☐	☐

- Switch language from Japanese to international (English-US)

ACTION	MRC	MRSC
Select Options -> Configuration in the main menu of the Meduser platform	☐	☐
Double-click the Regional and Language Settings icon	☐	☐

ACTION	MRC	MRSC
On the first tab card (Regional Options) , change the settings as shown in the picture below (US English).	☐	☐
Fig. 17: 		

ACTION	MRC	MRSC
On the second tab card (Language), change the settings as shown in the picture below (US English). Fig. 18: 	☐	☐
First press OK, next Apply. Remark: Pressing "OK" in the "Regional and Language Options" tab card starts an automatic reboot. Wait for two steps with this action.	☐	☐
Open the Windows Explorer (it could be necessary to switch to "Advanced User" to open the Explorer)	☐	☐
Delete the file C:\MedCom\config\Mri\MRSystem\MrVCheck.ref	☐	☐
Press the OK button and accept all messages, which will be shown with "YES" or "OK." A reboot follows!	☐	☐

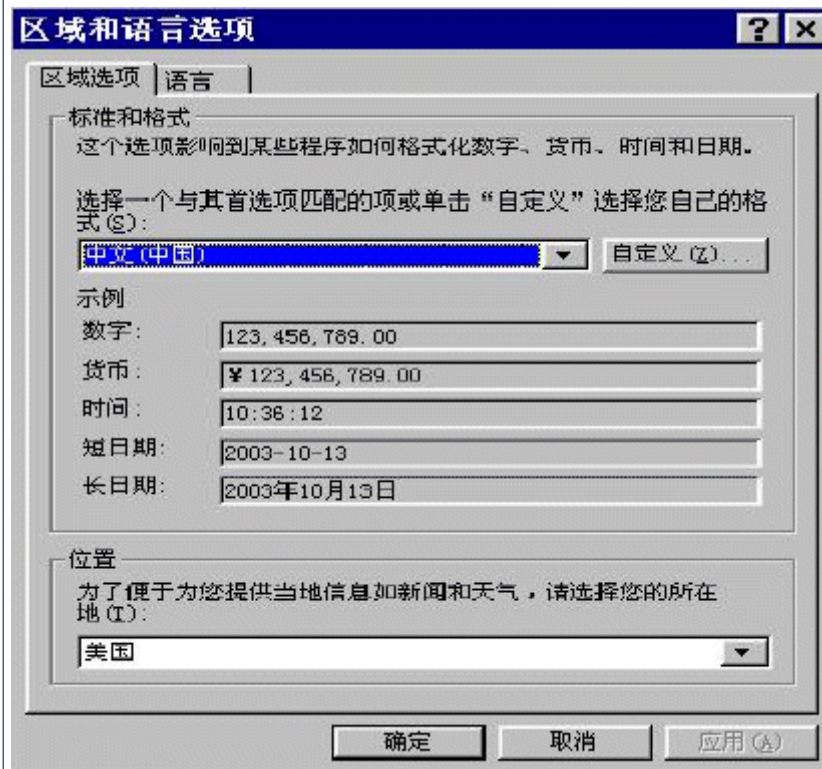
5.1.3 Language Settings Chinese

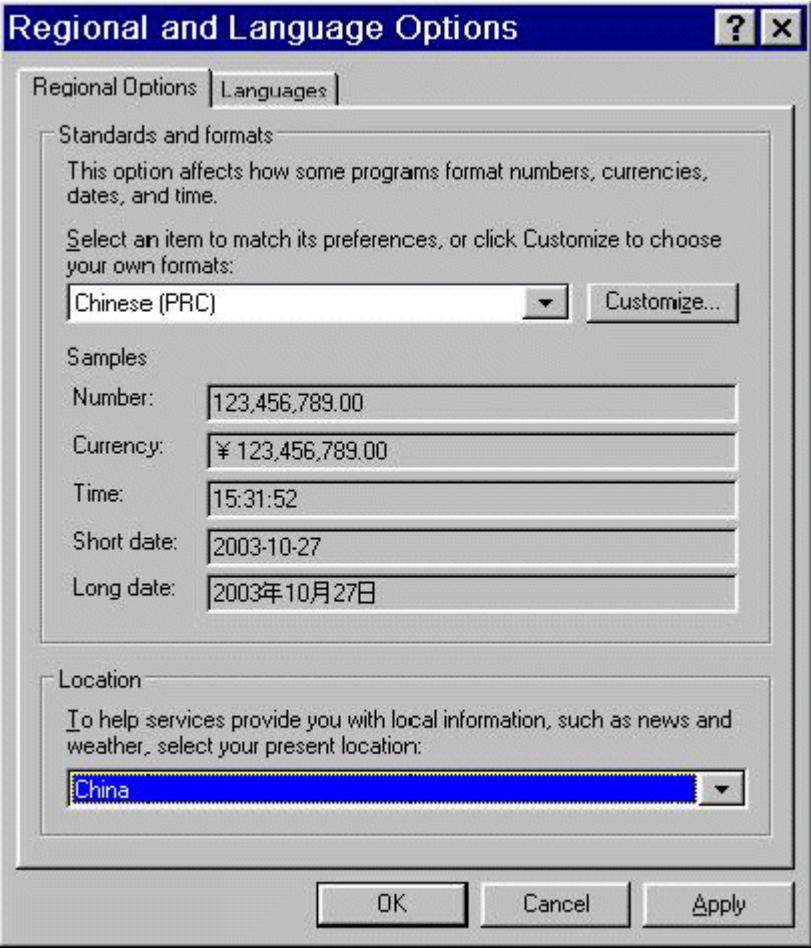
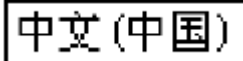
- Switch the language from international (US English) to Chinese

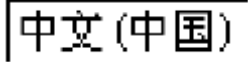
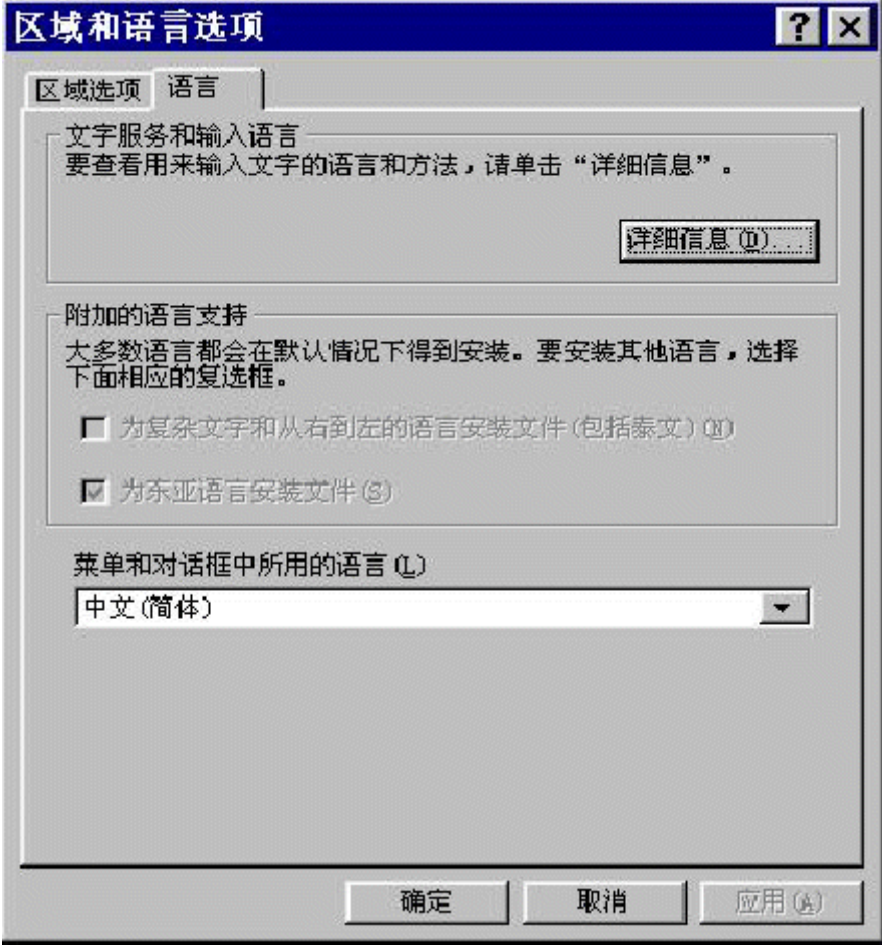
ACTION	MRC	MRSC
Select Options -> Configuration in the main menu of the Meduser platform	☐	☐

ACTION	MRC	MRSC
Double-click the Regional and Language Settings icon	☐	☐
In the highlighted field, choose Output pop-up in Chinese, if it was already switched during the Windows XP installation.	☐	☐

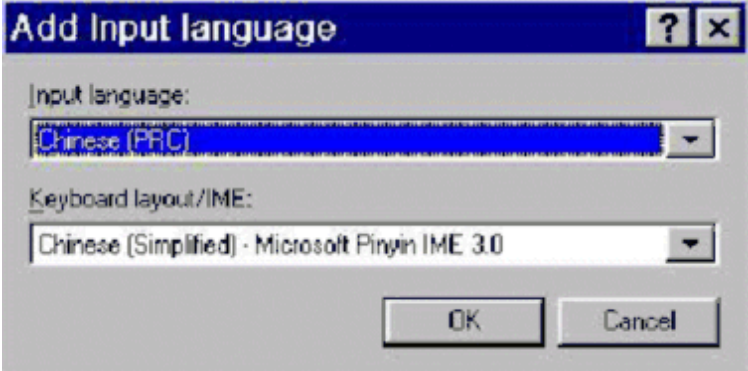
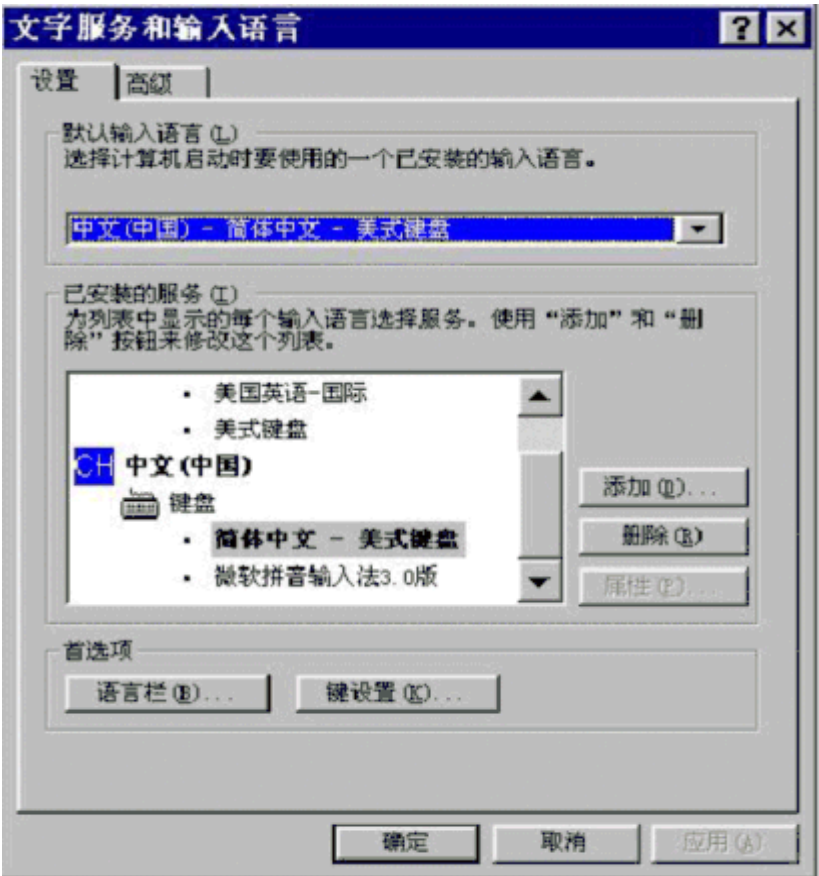
Fig. 19:

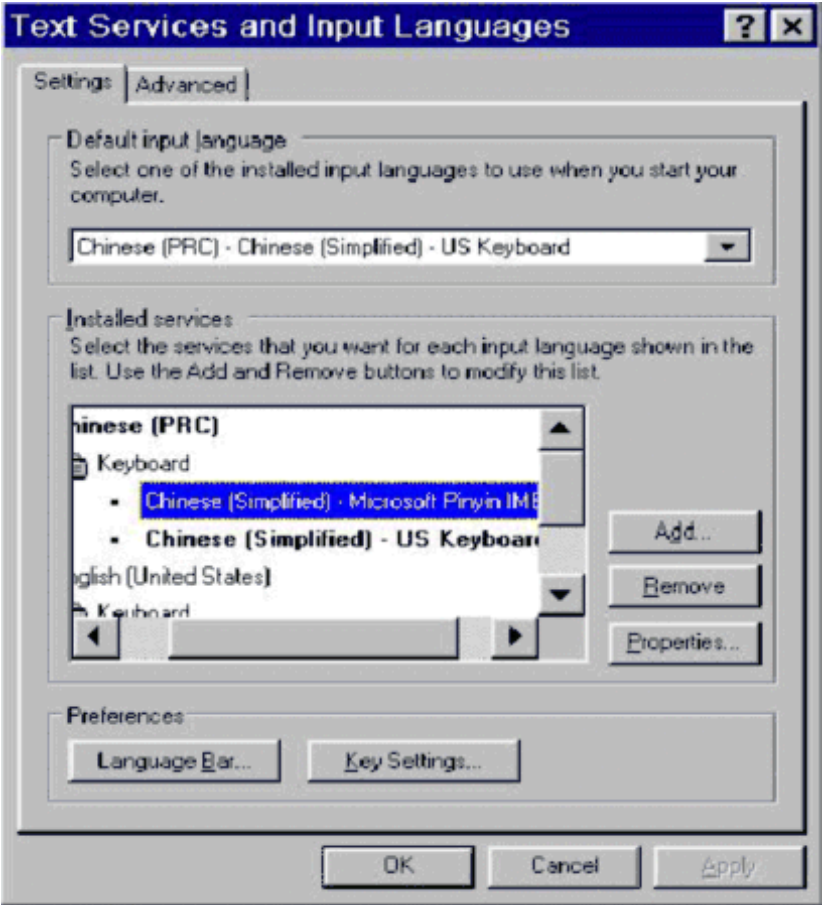


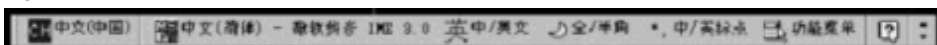
ACTION	MRC	MRSC
Otherwise, the international view of the pop-up is displayed. Select the entries shown below. Fig. 20:  The following item (see below) Fig. 21:  (Remark: scroll up, if the following symbol is not present)		

ACTION	MRC	MRSC
Switch to the 2. tab card "language" and choose the same Chinese symbol as in the first tab card Chinese symbol: Fig. 22:  Chinese view Fig. 23: 	☐	☐

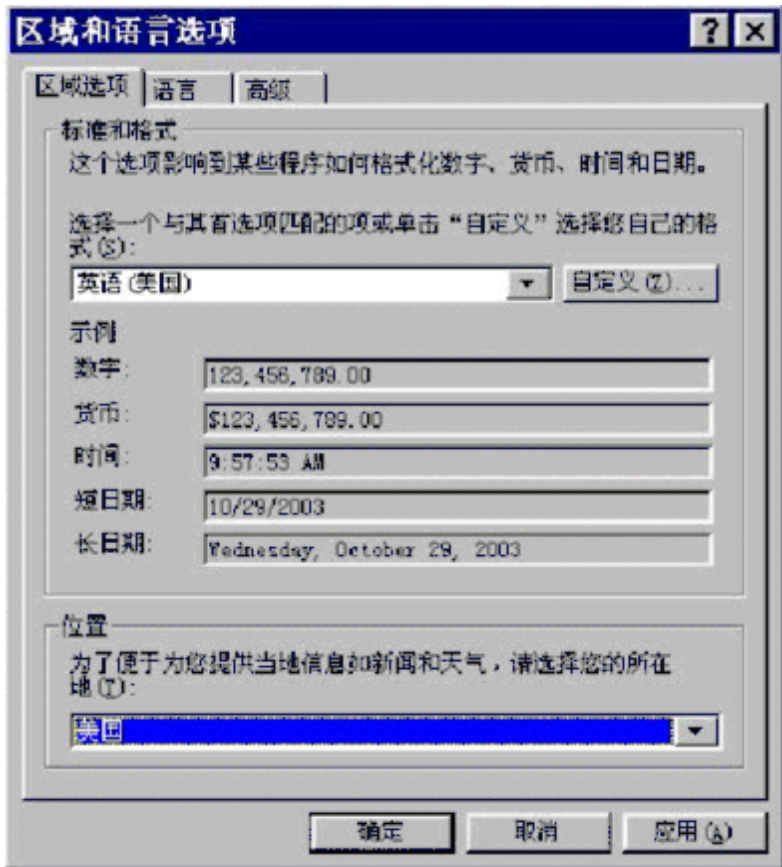
ACTION	MRC	MRSC
or international (US English) view Fig. 24:  The screenshot shows a Windows-style dialog box titled 'Regional and Language Options'. It has two tabs: 'Regional Options' and 'Languages'. The 'Languages' tab is active. Inside, there are three sections: 1. 'Text services and input languages' with a 'Details...' button. 2. 'Supplemental language support' with two checkboxes: 'Install files for complex script and right-to-left languages (including Thai)' (unchecked) and 'Install files for East Asian languages' (checked). 3. 'Language used in menus and dialogs' with a dropdown menu currently showing '中文(简体)'. At the bottom are 'OK', 'Cancel', and 'Apply' buttons.		

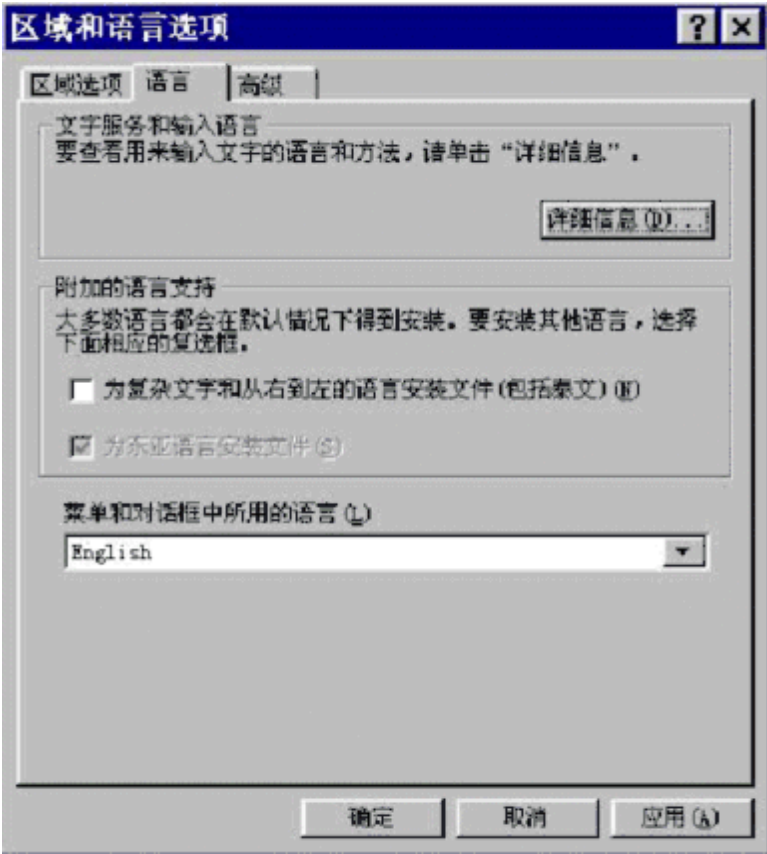
ACTION	MRC	MRSC
■ Click the "Details" button Note: If the necessary keyboard layout/IME is not present, press the "Add" button and select (see following pop-up) and press OK. Fig. 25: 	<input data-bbox="1305 275 1345 309" type="checkbox"/>	<input data-bbox="1409 275 1449 309" type="checkbox"/>
Subsequently, perform the following settings for your keyboard: Chinese view Fig. 26: 		

ACTION	MRC	MRSC
or international (US English) view Fig. 27: 		

ACTION	MRC	MRSC
Log in as meduser (if autologin is not enabled) after the reboot and wait until the main UI is up.	☐	☐
You can switch the output language of IME to Chinese by clicking it: (See the red symbol in the middle). Fig. 30:  or to English Fig. 31: 	☐	☐

- Switch the language from Chinese to international (US English)

ACTION	MRC	MRSC
Select Options -> Configuration in the main menu of the Meduser platform	☐	☐
Double-click the Regional and Language Settings icon	☐	☐
On the first tab card (Regional Options) , change the settings as shown in the picture below (US English). Fig. 32: 	☐	☐

ACTION	MRC	MRSC
On the second tab card (Language), change the settings as shown in the picture below (US English). Fig. 33: 	☐	☐
First press OK, next Apply. Remark: Pressing "OK" in the "Regional and Language Options" tab card starts an automatic reboot. Wait for two steps with this action.	☐	☐
Open the Windows Explorer (it could be necessary to switch to "Advanced User" to open the Explorer)	☐	☐
Delete the file C:\MedCom\config\Mri\MRSystem\MrVCheck.ref	☐	☐
Press the OK button and accept all messages, which will be shown with "YES" or "OK." A reboot follows!	☐	☐

5.2 Archiving

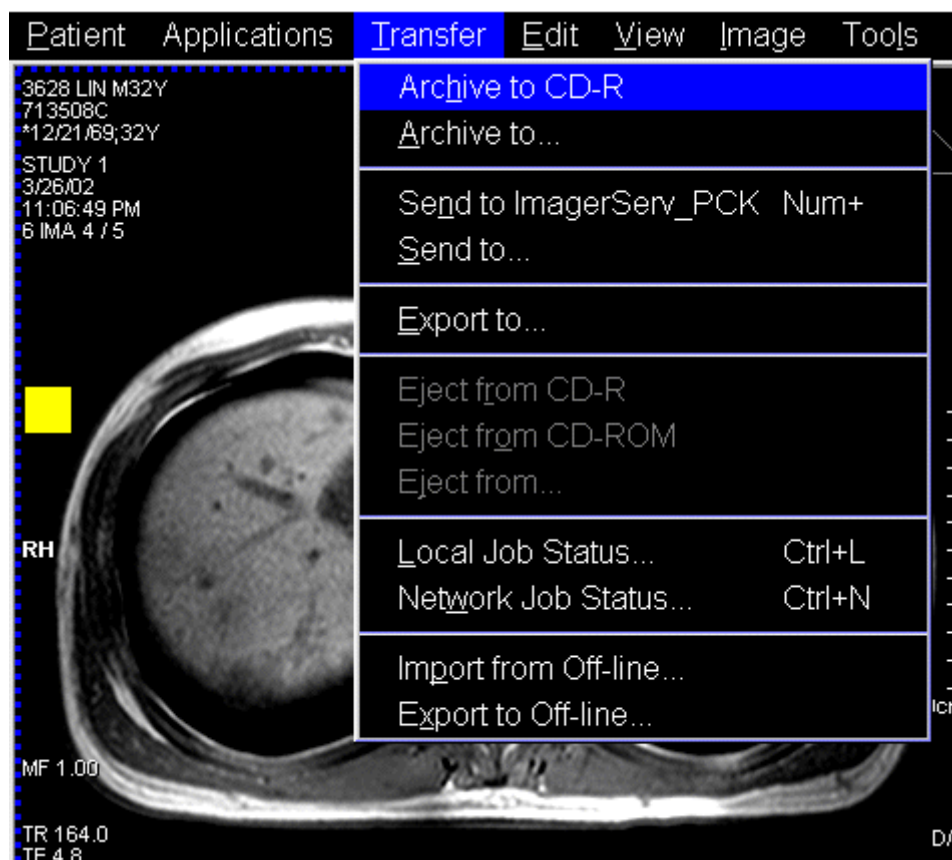
5.2.1 Configuring the CD-R drive as an Archive Device

If the customer requires the CD-R drive as an archive device, you can configure it as such in the service platform.

1. Start local service.
2. Enter the valid Service Key and confirm by clicking **OK**.
3. Select **Configuration**.
4. Under DICOM select **Offline Devices**.
5. Select the CD-R entry under the **Select logical name** pull-down menu.
6. Mark the checkbox **Archive device**.
7. Click **Save** and **Finish**.
8. Reboot by clicking **Home**.
9. Confirm the displayed windows by clicking **OK**.

After the system reboots, the following entry is displayed in the Viewer task card: "Transfer, Archive to CD-R".

Fig. 34:



5.2.2 Settings for Codonix Color printer NP1660M

5.2.2.1 Description:

Problems with the quality of the print out from Codonix Color printer NP1660M .

5.2.2.2 Workaround:

Check the settings of the printer if the configuration is set to color.

- Go in local service
- Select Configuration -> DICOM Print devices
- Go to step 2, HC device
- Select the entry for the Codonix Color printer under HC device
- Check if the following entry "color appearance" is set to **Color 24 Bit** if printed on paper (Medium type is set to paper).

Is Medium type is set to blue film, the entry for "color appearance" must be set to **Grayscale 8 bit** .

Is the setting is wrong, change it.

- Don't forget to save any changes made.

5.3 Configuration

5.3.1 Chorus

5.3.1.1 Prerequisite

Host and DICOM configuration for the protocol server is completed.

5.3.1.2 Configuration

Fig. 35: Chorus VD11A

Configuration | Reports | Home | Help

Chorus

Host Connection

Host: EH3XTX3S | 172.16.6.194
 Port: 9990

☐ Use secure connection (SSL)
☒ Use IP-Address in server call

Services

Connected Server:
☐ SLR RIS ☐ MLR RIS ☒ syngo.via

Protocol Repository Service: http://172.16.6.194:9990/ObjectCommunication/WEBSVCofsyngo.Re

☒ Enable Protocol Upload: http://172.16.6.194:9990/WebRisAdapter/Default.aspx
☒ Enable Protocol Download: http://172.16.6.194:9990/WebRisAdapter/Default.aspx
☒ Enable Protocol Admin: http://172.16.6.194:9990/WebRisAdapter/Default.aspx
☐ Enable automatic protocol update after software update

Options

☐ Ignore Modality type of downloaded protocol items in Exam Explorer path
☒ Use region for study description (e.g. region*examination)

PhoenixZIP

Generation

☒ Generate PhoenixZIP documents

Ready | reboot system pending | < Finish Save |

1. **Host:** List of available servers (for information only).
2. **Connected server:** select server type at your site.
3. **Protocol Repository Service:** URL of RIS web service.
4. **Protocol Upload:** URL to RIS parameter upload web UI.

5. **Protocol Download:** URL for distributed download and explicit download requests UI.
6. **Protocol Admin:** URL to RIS administrations UI.
7. **Software Update:** Enable/disable software update service.
8. **Ignore first path level in Dicom Location:** "distributed" protocols tree format to be displayed in "Exam Explorer".
9. **Web Application UIs on RIS side:** Chorus supports two technologies:
 - a) "Web Application" like the Novius implementation provides a web portal to display protocol requests from the MR modalities (**check box is set**).
 - b) Non web application like IWM; Windows Forms and Java Applications are automatically downloaded and automatically displayed (check box is empty) .
10. **Use region for study description:** Use region name for study description.

5.3.2 Chorus and PhoenixZIP

When **Chorus** is activated also **PhoenixZIP** will be activated (It's not possible to disable PhoenixZIP without deactivating Chorus).

5.3.2.1 Prerequisite

Host configuration for the protocol server is completed.

(Ask the RIS server administrator regarding URL's for the different services)

5.3.2.2 Configuration

Fig. 36: Chorus and PhoenixZIP

Configuration

Chorus

Host Connection

Host:

Port:

☐ Use secure connection (SSL)

☐ Use IP-Address in server call

Services

Protocol Repository Service:

☒ Enable Protocol Upload:

☒ Enable Protocol Download:

☒ Enable Protocol Admin:

☒ Enable automatic protocol update after software update

Options

☒ Ignore Modality type of downloaded protocol items in Exam Explorer path

☐ Use Web Page UI instead of Rich Client UI (e.g. Java Webstart)
Note: Siemens MLR RIS uses a Rich Client UI, Siemens SLR RIS uses a Web Page UI

☒ Use region for study description (e.g. region*examination)

PhoenixZIP

Generation

☒ Generate PhoenixZIP documents

Ready | < | Finish | Save |

Chorus**1. Host Connection**

- **Host:** Select host in the pulldown menu
- **Port:** Enter the port of the Chorus server
- **Use secure connection (SSL)** Set the checkbox if the server uses https
- **Use IP-Address in server call** Set if DNS is not used in the LAN

2. Services

- **Protocol Repository Service:** Enter the tail of the URL (When this URL is entered and saved the MR system will activate Chorus)
- **Enable Protocol Upload** Set checkbox if used and enter the tail of the URL
- **Enable Protocol Download** Set checkbox if used and enter the tail of the URL
- **Enable Protocol Admin** Set checkbox if used and enter the tail of the URL
- **Enable automatic protocol update after software updated** Set checkbox if allowed (Multi Language RIS is not supported by this software version)

3. Options

- Ignore Modality type of downloaded protocol items in Exam Explorer path (Eliminates the Prefix "MR@" when checkbox is set)
- Ignore Modality type of downloaded protocol items in Exam Explorer path
- Use region for study description (e.g. region examination)

PhoenixZIP

1. Generation

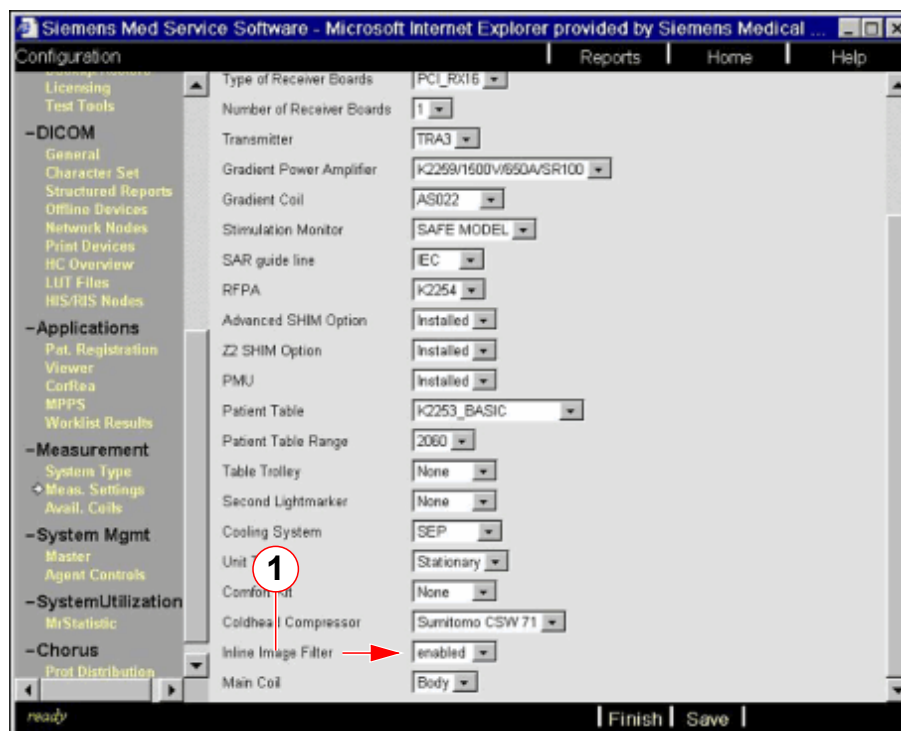
- **Generate PhoenixZIP documents** PhoenixZIP is automatically activated when Chorus is activated.

5.3.3 Enable/Disable Inline Image Filter

The following procedure describes how to do enable or disable the Inline Image Filter

- Open local service.
- Select Configuration.
- Under Measurement select Meas. Settings
- In the Measurement Settings Page go to the configuration entry Inline Image Filter and select enabled or disabled in the pull down menu.

Fig. 37: Service UI Measurement Settings



(1) Inline image filter

- After the changes have made, press **Save** and wait until **ready** and **reboot system pending** is displayed in the status line of the service UI.

- Perform the reboot by pressing **Finish** and **Home**.
Confirm any displayed reboot popup windows.

5.4 Troubleshooting

5.4.1 Customer network IP conflict with measurement network

If the customer local area network uses the same IP address range as the measurement network (host, MPCU, and Imager) the addresses for the measurement network must be changed, see following procedure:

- Start the local Service at the mrc, select **Configuration**
- Select **TCP/IP**
- Select the network adapter "[00000001] Measurement (internal)" this NIC is for used the for measurement network.
- Select the IP address, click the **Del** to remove the original entry.
- Add the new IP address (e.g. 192.168.4.1) by inserting it into the appropriate field and subnet mask (e.g. 255.255.255.248), click **Add** .
- Click **Save** to save the changes.
- Reboot the mrc
- Login as "meduser"
- Shutdown and switch off the system, wait ~20 seconds, boot the system again.

Now the measurement network should work with the new (measurement network addresses)

5.4.2 Problems during syngo Installation

5.4.2.1 Situation:

You are prompted to insert syngoMR B15 DVD 1/3 while the system is rebooting due to an interruption during installation. However, after inserting the CD, you are prompted to remove the disk because installation is supposedly already in progress.

5.4.2.2 Workaround:

1. Insert the CD labeled "syngoMR B15 DVD 1/3" and reboot the system.
2. When the syngo installation screen appears with the error message, press "CTRL + C" to cancel installation.
3. The DOS prompt appears on the left side of the screen. If the drive letter is "C:\", change to drive A:\ by entering "a:" and pressing <ENTER>.
4. Change the directory using "cd utils".
5. Type "fdisk" and press <ENTER> to start the partitioning program.
6. Select the 4th menu item to see a listing of existing partitions. Check the file system type of the first partition.
7. Press <ESC> to exit the current menu.
8. Select the 3rd menu item to delete partitions.

9. Depending on the partition type, select one of the items 1 through 4 to delete your first partition. If you select the wrong menu item, press "ESC" to return to step 8. This also applies if the message "NO ... PARTITION TO DELETE" appears after choosing a menu item.
10. After choosing the correct type for your first partition, the warning "DATA IN THE ... PARTITION WILL BE LOST" and the message "WHAT ... PARTITION DO YOU WANT TO DELETE [1]" appears.
Press <ENTER>. If you are prompted for the partition name, enter "MED_SYSTEM" (this will also be displayed on the top of your screen).
11. Press "Y" for **Yes** (depending on your keyboard, this could be the "Z" key).
12. The message "...PARTITION DELETED" appears.

Close fdisk by pressing <ESC> twice, then press a key to reboot.

The installation restarts, provided that your "syngoMR B15 DVD 1/3" is still inserted.

5.4.3 Workaround for Satellite Console connection problem

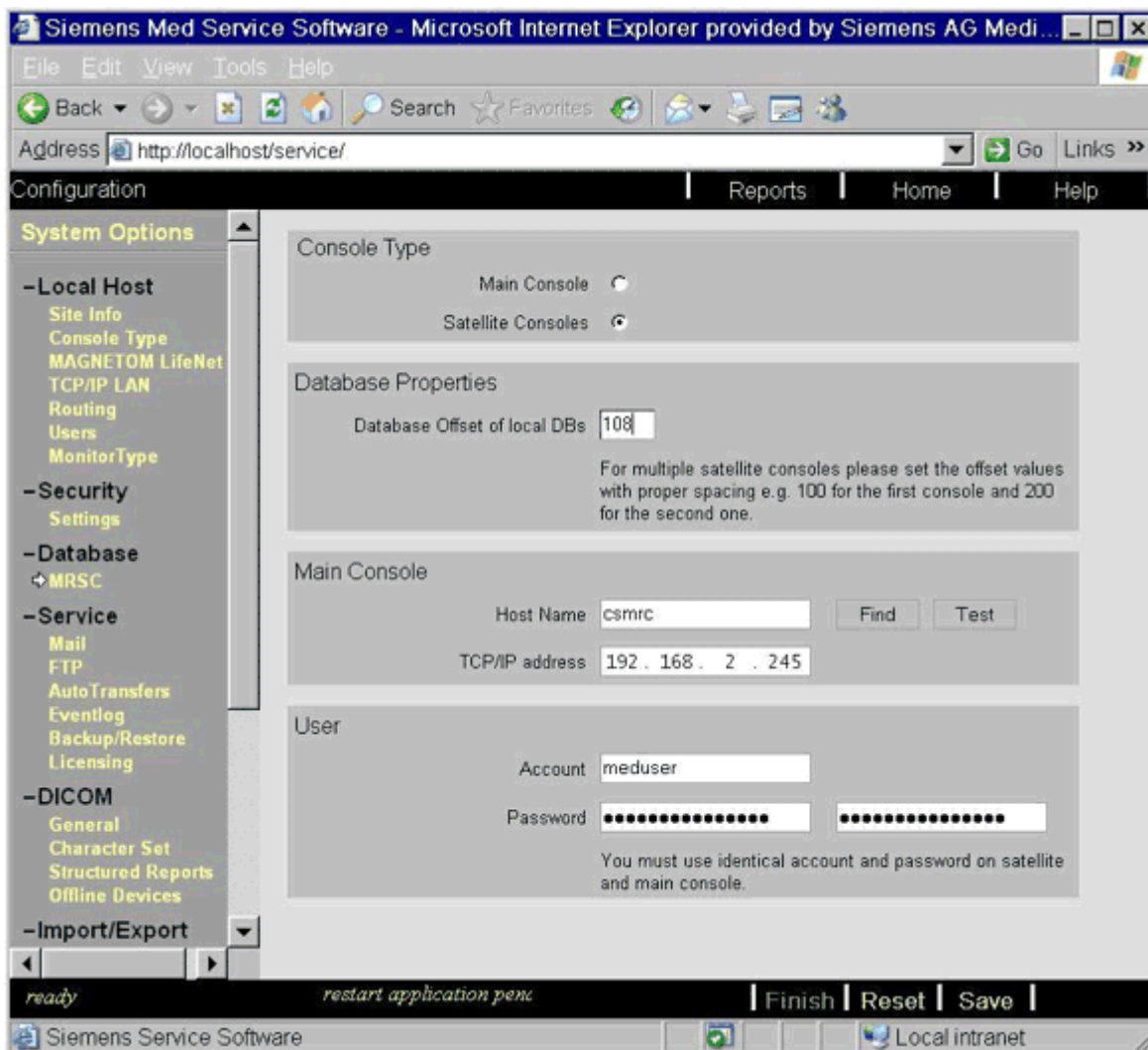
- Open local service
- Select Configuration
- Under "Database" **MRSC**:
 - Select **Remote Database** as "Database Type".
 - Under "Database Properties", select the value 102 for **Database Offset of local DB's** if only one satellite console is connected to the main console with a shared DB. If more than one satellite console is using the shared DB of the main console, each **Database Offset of local DB's** has to be assigned a different value, e.g., 102 for the first, 202 for the second., and so on.

Note: If the configuration of DBID is not unique within the network (configuration of satellite console), there is a risk that duplicate LOIDs will be created. This leads to problems in Filming and Archiving.

- Keep the other settings under Main Console and the entry "meduser" under the point "user" in the "Account" field.
- Enter meduser1 under password.

- Click **Save**.

Fig. 38:



- A reboot is initiated automatically.
Confirm the reboot popup window with "Yes"
The reboot starts (this will take about 5 min. until the reboot starts). In some cases a message popup "A component refused termination" appears. Confirm this popup with "OK".

5.4.4 Problems with MRSC Start up

5.4.4.1 Description:

The MRSC with shared database is not starting up. The message tab card appears with "No Access to MRC"

5.4.4.2 Workaround:

- On the MRC:
 - Go in local service
 - Select Configuration -> Database -> satellite console
 - Redo the configuration:

Select the entry for the satellite console in the pull-down menu "Select satellite console", delete this entry. If you are not able to delete it, go in "Service -> Ftp", select the satellite console entry under "Hosts" and delete the entry. Now go back to the configuration "Database -> satellite console" and add the host name and the IP address of the satellite console.

Use as password meduser1. Type this in under password and confirm.
 - Press Save
 - Continue with Finish.
 - Close local service
- On the MRSC
 - Go in local service
 - Select Configuration -> Database MRSC
 - Select **Remote Database** as "Database Type" .
 - Select the entry of the main console and delete it in the pull-down menu **Select main console** under "Main console" . If you are not able to delete it, go in "Service -> Ftp", select the main console entry under "Hosts" and delete the entry. Go back to "Database MRSC"
 - Under "Database Properties" **Database Offset of local DB's**, use the value 100 if only one satellite console is connected to the main console with shared DB. If more than one Satellite console are using the shared DB of the main console, each **Database Offset of local DB's** must have a different value e.g. 100 for the first, 200 for the second ...
 - Add the entry of the main console (host name and IP address) in the pull-down menu **Select main console** under "Main console" .
 - Use meduser1 as password under Users password . Type this in under password and confirm.
 - Press **Save**
 - Reboot manually if no automatic reboot is performed.

Chapter Site-specific data: Information added how to Backup/Restore site specific data on external USB device.

There are no Hazard IDs in this document.

All documents may only be used for rendering services on Siemens Healthcare Products. Any document in electronic form may be printed once. Copy and distribution of electronic documents and hardcopies is prohibited. Offenders will be liable for damages. All other rights are reserved.

healthcare.siemens.com/services

Siemens Healthineers Headquarters
Siemens Healthcare GmbH
Henkestr. 127
91052 Erlangen
Germany
Telephone: +49 9131 84-0
siemens.com/healthineers

Print No.: MR-000.816.24.02.02 | Replaces: MR-000.816.24.01.02
Doc. Gen. Date: 07.08 | Language: English
Part No.: n.a.
© Siemens Healthcare GmbH, 2007

siemens.com/healthineers