

MAGNETOM

Software System

Software Installation syngo MR B19 (VB19B)

Document Version

Siemens reserves the right to change its products and services at any time.

In addition, manuals are subject to change without notice. The hardcopy documents correspond to the version at the time of system delivery and/or printout. Versions to hardcopy documentation are not automatically distributed.

Please contact your local Siemens office to order a current version or refer to our website <http://www.healthcare.siemens.com>.

Disclaimer

Siemens provides this documentation “as is” without the assumption of any liability under any theory of law.

The content described herein requires superior understanding of our equipment and may only be performed by qualified personnel who are specially trained for such installation and/or service.

Copyright

“© Siemens, 2017” refers to the copyright of a Siemens entity such as:

Siemens Healthcare GmbH - Germany

Siemens Aktiengesellschaft - Germany

Siemens Shenzhen Magnetic Resonance Ltd. - China

Siemens Shanghai Medical Equipment Ltd. - China

Siemens Healthcare Private Ltd. - India

Siemens Medical Solutions USA Inc. - USA

Siemens Healthcare Diagnostics Inc. - USA and/or

Siemens Healthcare Diagnostics Products GmbH - Germany

1	General	5
1.1	syngo MR system software	5
1.1.1	Introduction	5
1.1.2	Prerequisites	5
1.1.3	Abbreviations	5
2	Installing Software (MRAWP(MRC) / MRWP(MRSC))	6
2.1	Prerequisites	7
2.2	BIOS Update	8
2.2.1	Load BIOS defaults for HP Z420	8
2.2.2	Load BIOS defaults for HP Z440	8
2.3	Install Software from system recovery image	10
2.4	Final checks	12
3	Site-specific Data	13
3.1	Site-specific data	13
3.1.1	Backup site-specific data on CD	13
3.1.2	Backing up site-specific data on external USB device	14
3.1.3	Create COBR Image Backup	14
3.1.4	Restore of the COBR Image	37
3.1.5	Restore site-specific data from CD	78
3.1.6	Restore site-specific from external USB device	79
4	Configuration	82
4.1	Configuration	82
4.1.1	General	82
4.1.2	Network Configuration Restrictions	83
4.1.3	Activation of Security Settings	83
5	Additional information, troubleshooting	84
5.1	Appendix	84
5.1.1	Monitor Settings	84
5.1.2	License update	85
5.1.3	Restoring Customer Protocols from Last Backup	87
5.1.4	Dongle exchange and license update	87
5.1.5	Customer protocols	89
5.1.6	Language settings and font size	89
5.1.7	Keyboard and language setting	90
5.1.8	Enable/Disable Inline Image Filter	127
5.1.9	Reset MRVCheck	127
5.1.10	File cleanup for image creation	128
5.2	Hints	129
5.2.1	Exporting and backing up customer measurement programs	129
5.2.2	Restoring and importing customer measurement programs	129
5.2.3	Information for the account administrator	129
5.2.4	Changing network addresses for the measurement network	130

5.2.5	Use of customer-built coils or 3rd party coils with Magnetom systems	130
5.2.6	Startup error "Found Databases with Equal IDs" on MRSC	131
5.2.7	Drive letter assignment for the CD/DVD ROM drive and/or for the CD/DVD-R drive	131
5.2.8	Drive letter assignment for hard disks.	134
5.2.9	COBR Backup was interrupted, no System Boot	137
6	Changes	138
7	List of Hazard IDs	139

1.1 syngo MR system software

1.1.1 Introduction

The installation and service of equipment described herein is to be performed by qualified personnel employed or authorized by Siemens Medical Solutions or one of its affiliates. Assemblers and other personnel not employed by nor directly affiliated with nor authorized by Siemens Medical Solution or one of its affiliates` technical services are directed to contact the local office of Siemens Medical Solutions or its affiliates before attempting installation or service procedures.

1.1.2 Prerequisites

To configure the syngo-based MR system, the following training is required in combination with this software manual:

- syngo training / syngo MR training
- Network training (TCP/IP)
- DICOM training



All passwords required to install and configure the system are available in the "CS Knowledge Base" <https://smskb.medical.siemens.com/selfservice> under "Password List" (except for the service key (variable password), which is handled by the districts).

1.1.3 Abbreviations

EDB	Event Data Base
n.a.	Not applicable
SRS	Siemens Remote Service
QDB	Quality Assurance Data Base

This chapter describes how to install software with the COBR tool.



Each host (MRC or MRSC) has its own individual set of recovery image DVD's included with the system delivery.



No second monitor is supported in VB19B. Please be sure that only one monitor is connected to MRAWP or MRWP.

2.1 Prerequisites

Required media:

- Host (MRC or MRSC) specific recovery image DVD
- BIOS CD

2.2 BIOS Update

The BIOS update optimizes the setting for NUMARIS/4 for your host type, sets the NUMA-RIS/4 BIOS password and boot device order

2.2.1 Load BIOS defaults for HP Z420

- Reboot the system (power on the Host).
- During the boot, press **F10** to access the BIOS setup.
- Enter the BIOS-Password.
- Check date / time, correct if necessary.
- Select **Storage**
- Highlight the **Boot Order** and press **Enter**.
- Select **Optical Drive** press **Enter** to activate "Drag mode".
- Move the "Optical Drive" entry to the top of the list.
- Press **Enter** to "Drop" the setting.
- Press **F10** and accept "Yes" by pressing **Enter**.
- Select **File**, select **Save Changes and Exit**
- Press **F10** and accept "Yes" by pressing **Enter**.
- The BIOS-update will run automatically.
- Remove the BIOS-CD after the automatic eject.
- Restart the system by pressing **Ctrl + Alt + Del** simultaneously.

2.2.2 Load BIOS defaults for HP Z440

- Reboot the system (power on the Host).
- During the boot, press **F10** to access the BIOS setup.
- Enter the BIOS-Password.
- Check date / time, correct if necessary.
- Select **Advanced**
- Select **Boot Options** .
- Move the scroll bar to the bottom
- Select **Legacy Boot Order** and press **Enter**.
- Select the entry **CDROM** and press **Enter** the to activate "Drag mode"
- Move the **CDROM** entry to the top of the list.
- Press **Enter** to "Drop" the setting.
- Press **Esc** to leave the menu
- Press **Esc** to leave the menu
- Select **Main**, select **Save Changes and Exit**
- Select **Save Changes and Exit** and press **Enter**.

- Select "Yes" by pressing **Enter**.
- The BIOS-update will run automatically.
- Remove the BIOS-CD after the automatic eject.
- Restart the system by pressing **Ctrl + Alt + Del** simultaneously.

2.3 Install Software from system recovery image



All data will be lost!

- Power off the system on the alarm box.
- Disconnect the power line on the mars.
- Switch on the system on the alarmbox.
- During the boot, press **F10** to access the BIOS setup.
- Enter the BIOS-Password.
- Check date / time, correct if necessary.
- Select **Advanced**
- Select **Boot Options** .
- Move the scroll bar to the bottom
- Select **Legacy Boot Order** and press **Enter**.
- Select the entry **CDROM** and press **Enter** the to activate "Drag mode"
- Move the **CDROM** entry to the top of the list.
- Press **Enter** to "Drop" the setting.
- Press **Esc** to leave the menu
- Insert the COBR Backup Restore DVD in the DVD-ROM drive
- Select **Main**, select **Save Changes and Exit**
- Select **Save Changes and Exit** and press **Enter**.
- Select "Yes" by pressing **Enter**.
- The COBR tool will boot (several cmd windows will be displayed, do not interact).
- In the "Backup/Restore UI" Select **Restore System**, click **OK**
- The system recovery will proceed automatically,
- Eject the drive when requested; insert the next DVD and click **OK**
- After completed restore from the last media answer the question "**If you want to insert the next DVD for Restore please press the 'Yes' button, otherwise press 'NO'**" with **NO**
- The partitions will now be restored from the image files.
- Remove the Recovery image DVD, click **OK**, re-boot will follow automatically
- Now the COBR post restore process is started. Note:



Please don't klick or hit any key ! Just wait.

- Wait until the pop up "**Preparation for mars installation finished. Please start the MARS and press OK**" is displayed.
- Connect the power line on the mars and confirm the popup with **OK**.
- Wait until the MARS Installation has finished and the syngo main UI is up and running.

- Perform a reboot of the host.
- Insert BIOS-CD Vxxxx xxxxx-MRxxxBIOS-CD in to the DVD/CD-ROM drive.
- The BIOS-update will run automatically.
- Remove the BIOS-CD after the automatic eject.
- Restart the system by pressing **Ctrl + Alt + Del** simultaneously.



"Autologin" is enabled by default in the COBR image.

2.4 Final checks



Check if hotfixes and Service Packs are available for this software version according to various Speedinfos and Update Instructions.

- Restore the most recent "site specific data" (e.g. customer protocols).
- Perform the "Quality Assurance" procedure: "Image Orientation" in the service platform.
- Make sure that the system is producing images.
- Check that the system time, time zone and date is correct.



Please contact the customer / admin if the advanced user should be enabled; also ask the customer / admin what shall be set as advanced user password and for the other users.

3.1 Site-specific data

3.1.1 Backup site-specific data on CD

Required time: 10-60 min

Required tools / parts: CD-R or USB device



Use only CD-R 's with medical grade, e.g. Siemens Medical CD-R 80 or TDK Medical CD-R.

Save the site-specific data each time you make changes to the site-specific configuration, tune-up, or sequences and protocols. In case of a hard disk crash, you will be able to restore a valid set of site-specific data.

- Login as **administrator**
- Start the service platform, select **Backup & Restore**.
- Insert a recordable CD-disk in the CD/DVD-burner.
- Select **Backup** in the "Command" pull-down menu.
- In the menu Drives, select the target drive CD/DVD-burner (drive letter **S:**) used to save the data.
- Select **Master-Backup** in the "Packages" menu.
- Click **Go** to start the "Backup", check for errors in the log-window after the backup has been completed.

Setting that are not part of the site-specific data:

Remote access mode (No Access, Limited Access, etc. in the customer UI).

Auto-transfer flag ("Configuration", "Service", "Auto Transfers")

Postscript Printer Configuration (Customer paper printer)

ImageText Settings (ImageText Settings from ImageText Dialog)

Time zone

Transfer Rules

Autologin for meduser (frontend username)

Frontend username (e.g. **meduser**)

Secure DICOM printer setting (must be manually configured)



During backup at the MRSC, a "**CorRea.bru**" error may be displayed, please ignore this error.

3.1.2 Backing up site-specific data on external USB device



Please check the external media (e.g. USB-disk) first with the virus scanner of the service PC. Be sure that the external media used is free of viruses.

- Login as **administrator** ↵.
- Start the "Internet Explorer"
- Connect the USB device.
- Check what drive letter was set for the USB device with the windows file explorer.
- Login to the service platform, select **Backup & Restore**.
- Select Backup in the Command menu, select drive letter of the USB device e.g. **F** in the Drives menu..
- Select **Master-Backup** in the Package menu, start the backup by clicking **GO**.
- (Estimated time ~10-30 min) check for errors using the log window.
- Click **Home** to leave Backup & Restore.

Setting that are not part of the site-specific data:

Remote access mode (No Access, Limited Access, etc. in the customer UI).

Auto-transfer flag ("Configuration", "Service", "Auto Transfers")

Postscript Printer Configuration (Customer paper printer)

ImageText Settings (ImageText Settings from ImageText Dialog)

Time zone

Transfer Rules

Autologin for meduser (frontend username)

Frontend username (e.g. **meduser**)

Secure DICOM printer setting (must be manually configured)



During backup at the MRSC, a "**CorRea.bru**" error may be displayed, please ignore this error.

3.1.3 Create COBR Image Backup

3.1.3.1 Manual Backup of Partition C Image to USB-Disk

This procedure describes how to back up the complete software system for MR AWP (MRC) / MR WP (MRSC) (all hard disks and partitions) to a USB-disk.





If you get an error popup with the message "Backup/Restore failed. The image*.tip is inconsistent" during the image consistency check from the COBR tool, please connect the USB device to the USB 2.0 port. With some USB devices it could be possible that you get this error message when the device is used on an USB 3.0 port. Repeat the image consistency check or repeat the image creation.



Prerequisites



- USB-disk with sufficient free space left.
- Autologin is disabled
- Password for syadmin is set and known
- Masterbackup is created and stored to partition ASIS drive letter F:.

Procedure

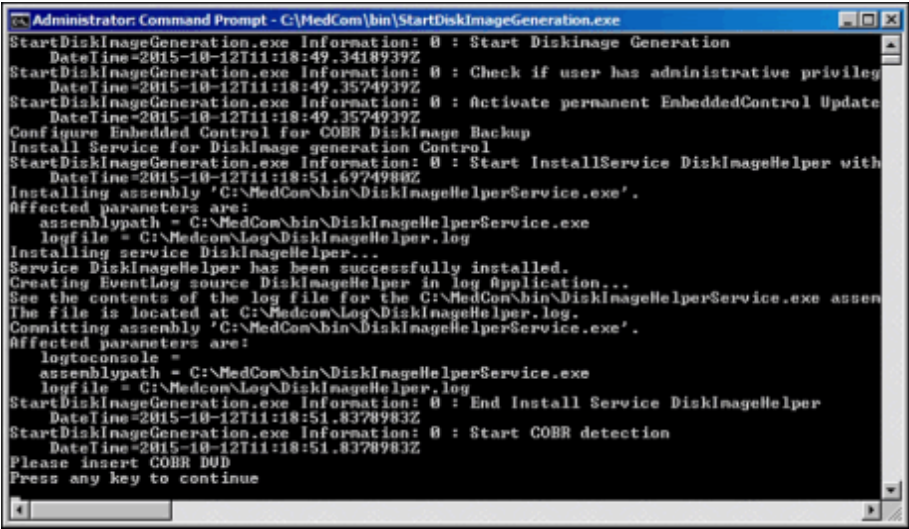ACTION	MRAWP/MRWP
■ Reboot and log in as syadmin	☐
■ Invoke a command shell (Start ⇒ All Programs ⇒ Accessories ⇒ Command Prompt)	☐
■ In the command shell, type in: <code>C:\MedCom\bin\StartDiskImageGeneration.exe</code> ↵	☐
<i>Fig. 1: Start COBR Image Creation</i> 	
The command above will switch EC to permanent update mode prepare the system for COBR image creation.	
■ The COBR tool DVD is requested. Insert the first DVD from the last created backup image into the CD/DVD-ROM drive.	☐
■ Press any key e.g. <SPACE> to start the COBR image backup.	☐

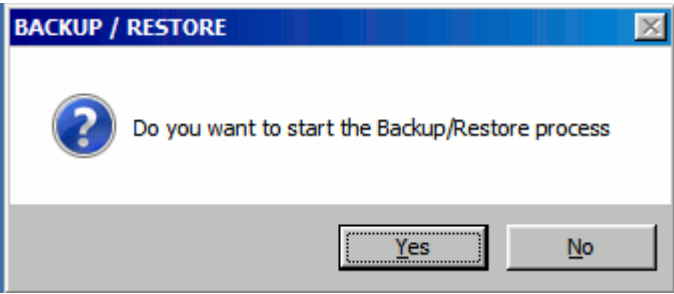ACTION	MRAWP/MRWP
Confirm the displayed popup with Yes. Fig. 2: Backup/Restore popup 	☐
The running script performs a reboot at the end. This may take some time. Please wait!.	☐
A reboot is performed and the COBR tool boots.	☐
Wait until the COBR tool window is displayed. Fig. 3: COBR main Backup/Restore popup 	☐
Connect the USB medium.	☐

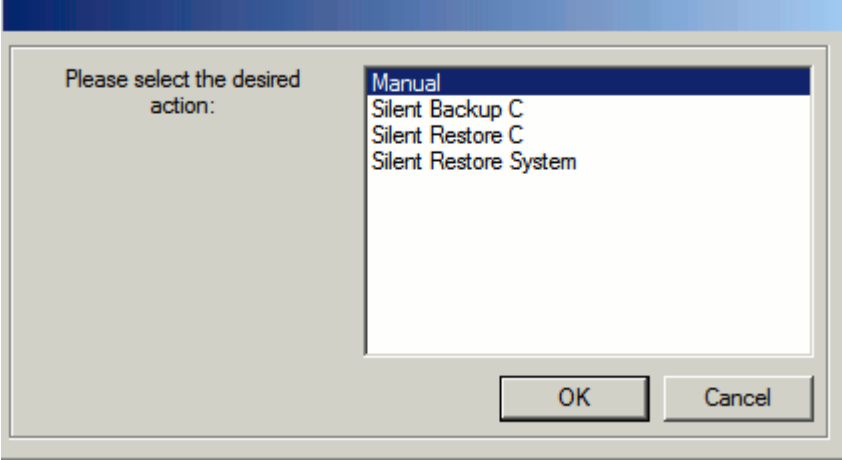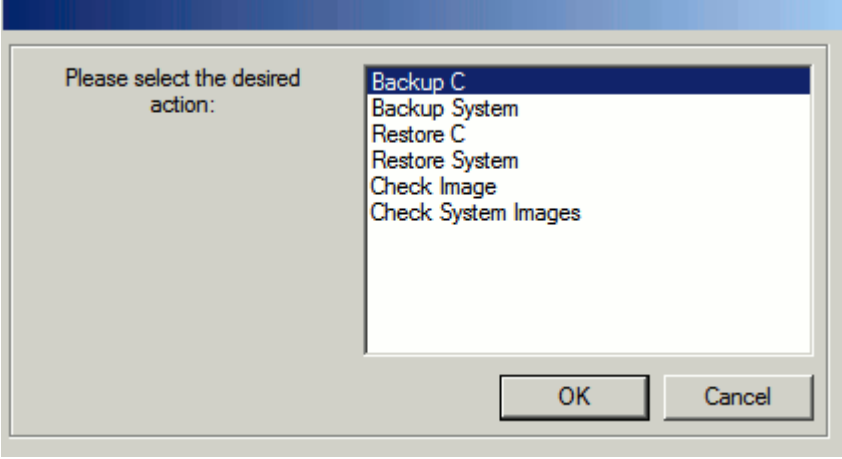
ACTION	MRAWP/MRWP
In the Backup/Restore popup select Manual and press OK Fig. 4: Manual Backup/Restore popup 	☐
In the next Backup/Restore popup select Backup C and press OK. Fig. 5: Backup/Restore popup - Backup C 	☐

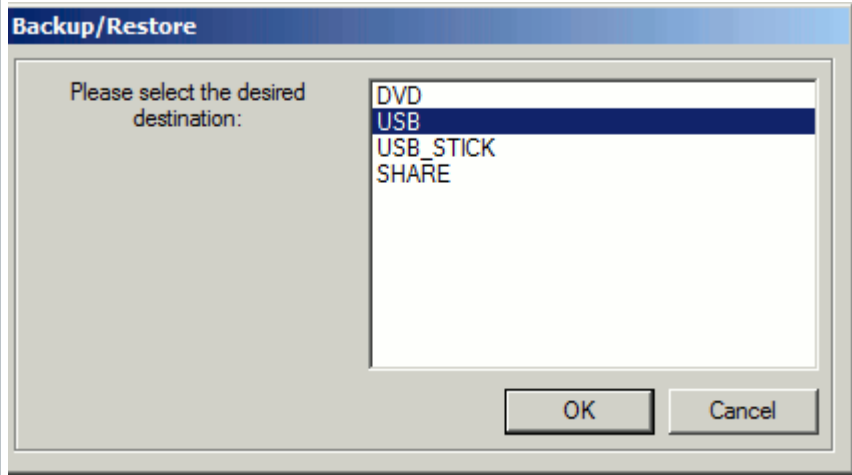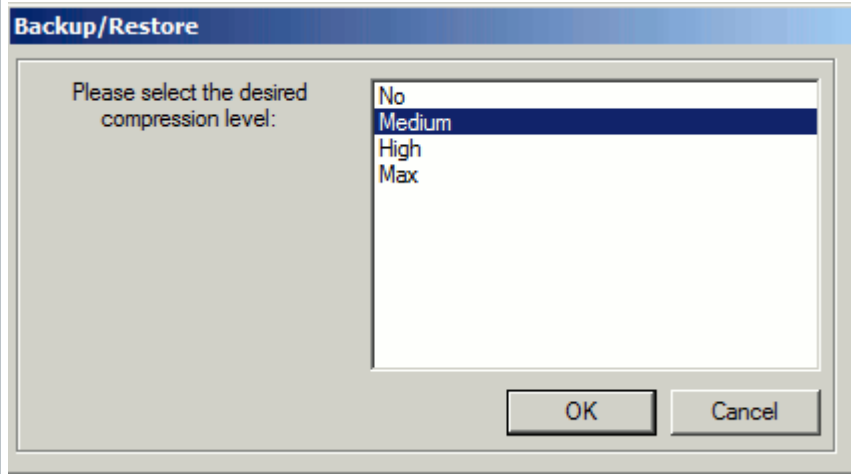ACTION	MRAWP/MRWP
In the next Backup/Restore popup select <code>USB</code> as the medium for the backup image and press OK. Fig. 6: Select medium USB 	☐
In the next Backup/Restore pop-up select the compression level <code>medium</code> of the data and press OK. Fig. 7: Backup/Restore popup - compression level 	☐

ACTION

- Now the system back up will run. The following is an example of the screen display.

Fig. 8: COBR Acronis image backup example output

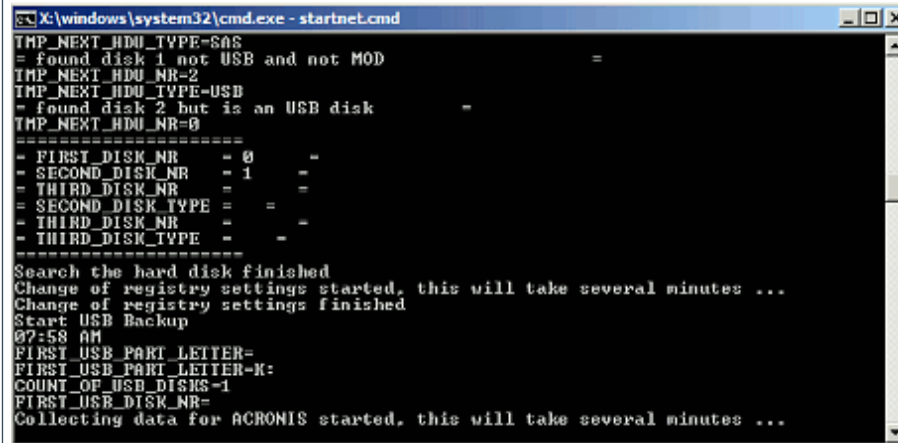

```

X:\windows\system32\cmd.exe - startnet.cmd
IMP_NEXT_HDD_TYPE=SAS
= found disk 1 not USB and not MOD
IMP_NEXT_HDD_NR=2
IMP_NEXT_HDD_TYPE=USB
= found disk 2 but is an USB disk
IMP_NEXT_HDD_NR=0
=====
- FIRST_DISK_NR = 0
- SECOND_DISK_NR = 1
- THIRD_DISK_NR =
- SECOND_DISK_TYPE =
- THIRD_DISK_NR =
- THIRD_DISK_TYPE =
=====
Search the hard disk finished
Change of registry settings started, this will take several minutes ...
Change of registry settings finished
Start USB Backup
07:58 AM
FIRST_USB_PART_LETTER=
FIRST_USB_PART_LETTER=M:
COUNT_OF_USB_DISKS=1
FIRST_USB_DISK_NR=
Collecting data for ACRONIS started, this will take several minutes ...


```

MRAWP/MRWP




- When the backups of all disks are written answer the popup "Do you want to create a bootable DVD" with "No".

Fig. 9: Create boot DVD pop-up

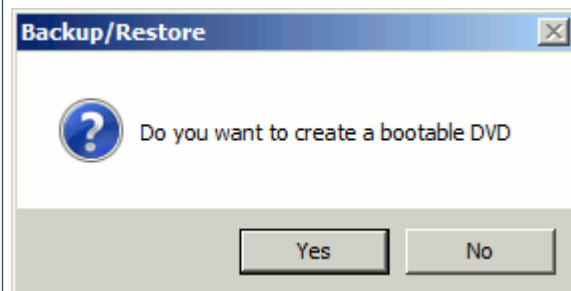

- Now a checksum (MD5) is created for each image file. This will take some minutes.

Fig. 10: MD5 checksum creation

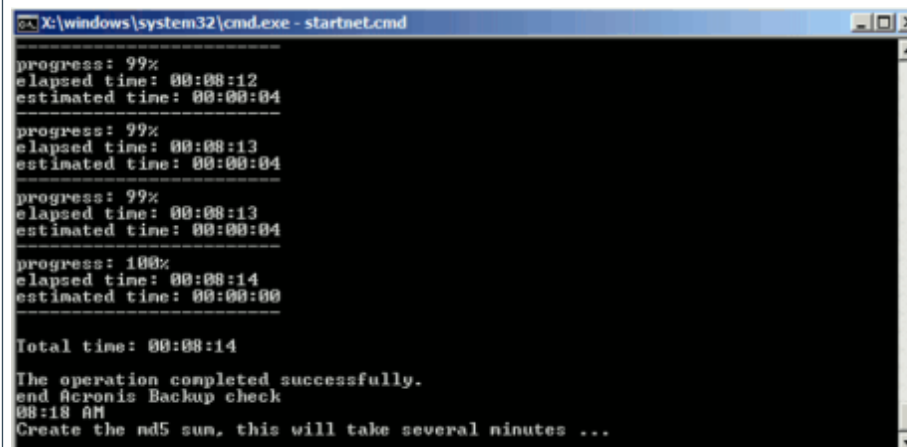





```

X:\windows\system32\cmd.exe - startnet.cmd
progress: 99%
elapsed time: 00:08:12
estimated time: 00:00:04
-----
progress: 99%
elapsed time: 00:08:13
estimated time: 00:00:04
-----
progress: 99%
elapsed time: 00:08:13
estimated time: 00:00:04
-----
progress: 100%
elapsed time: 00:08:14
estimated time: 00:00:00
-----
Total time: 00:08:14
The operation completed successfully.
end Acronis Backup check
08:18 AM
Create the md5 sum, this will take several minutes ...

```

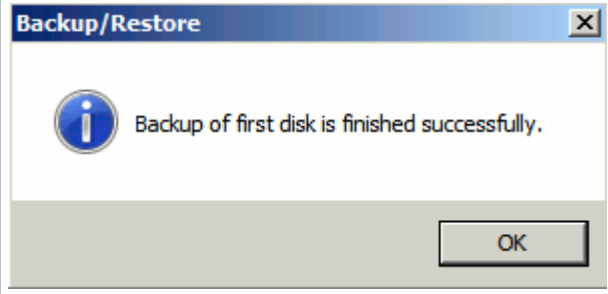

ACTION	MRAWP/MRWP
Confirm the success popup by pressing OK . Fig. 11: Backup finished successfully 	☐
The main Backup/Restore popup will appear. Fig. 12: COBR main Backup/Restore popup 	☐

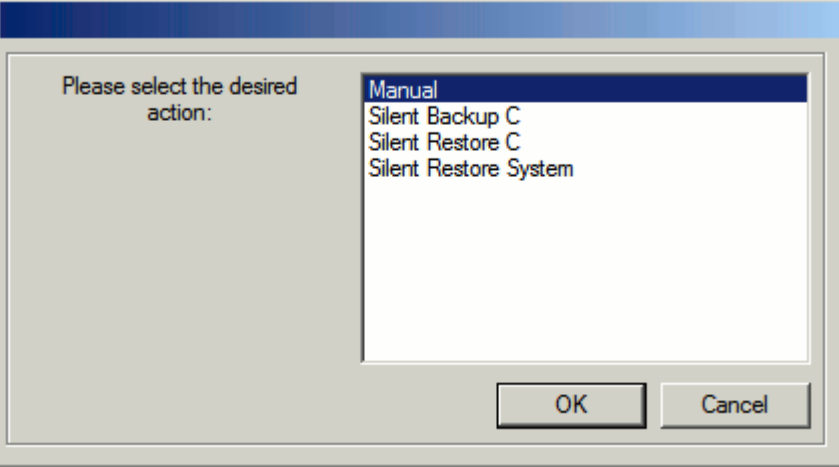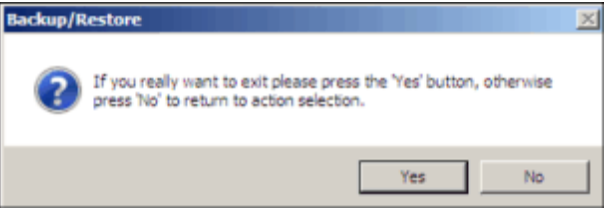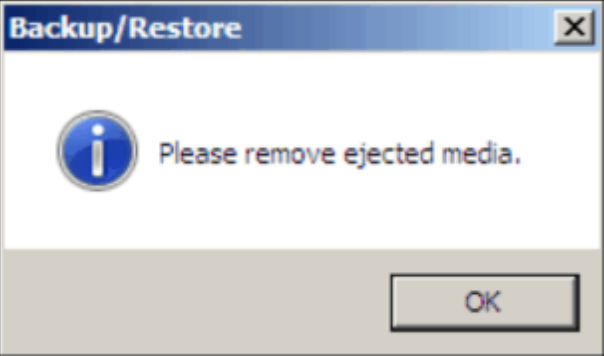ACTION	MRAWP/MRWP
Press Cancel to exit the Backup/Restore procedure. <i>Fig. 13: Main Backup/Restore popup</i> 	☐
Confirm the following popup with Yes. <i>Fig. 14: Confirmation popup</i> 	☐
The CD/DVD-ROM drive opens and the following popup is displayed. <i>Fig. 15: Backup/Restore remove ejected medium popup</i> 	☐
Remove the ejected disk from the drive.	☐
Close the drive and click OK to confirm.	☐
The system is rebooted.	☐
Disconnect the USB medium.	☐

ACTION	MRAWP/MRWP
After the reboot the windows login ist dispalyed. Please wait until a further reboot is performed (this takes about 30 seconds).	☐
The windows login is now displayed for the second time after the reboot. Now login as user meduser and wait until the main UI is displayed.	☐

3.1.3.2 Manual Backup of System Images to USB-Disk

This procedure describes how to back up the complete software system for MR AWP (MRC) / MR WP (MRSC) (all hard disks and partitions) to a USB-disk.





If you get an error popup with the message "Backup/Retore failed. The image*.tip is inconsistent" during the image consistency check form the COBR tool, please connect the USB device to the USB 2.0 port. With some USB devices it could be possible that you ge this error message when the decice is used on an USB 3.0 port. Repeat than the image consistency check or repeat the image creation.


Prerequisites



- USB-disk with sufficient free space left.
- Autologin is disabled
- Password for syadmin is set and known
- Masterbackup is created and stored to partition ASIS drive letter **F:**.




Procedure



ACTION	MRAWP/MRWP
Reboot and log in as syadmin	☐
Invoke a command shell (Start ⇒ All Programs ⇒ Accessories ⇒ Command Prompt)	☐

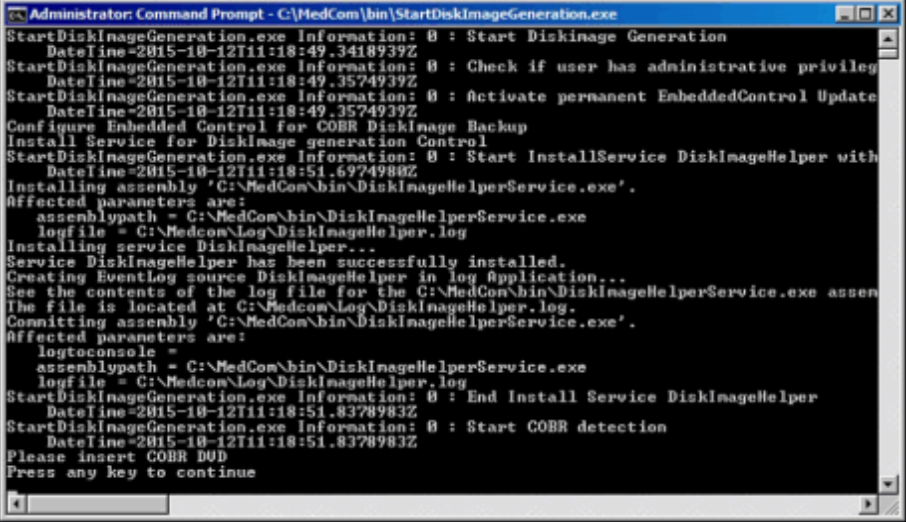
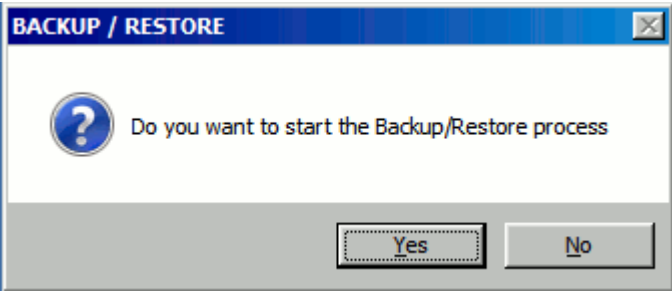
ACTION	MRAWP/MRWP
In the command shell, type in: <code>C:\MedCom\bin\StartDiskImageGeneration.exe ↵</code> Fig. 16: Start COBR Image Creation  The command above will switch EC to permanent update mode prepare the system for COBR image creation.	☐
The COBR tool DVD is requested. Insert the first DVD from the last created backup image into the CD/DVD-ROM drive.	☐
Press any key e.g. <SPACE> to start the COBR image backup.	☐
Confirm the displayed popup with Yes. Fig. 17: Backup/Restore popup 	☐
The running script performs a reboot at the end. This may take some time. Please wait!.	☐
A reboot is performed and the COBR tool boots.	☐

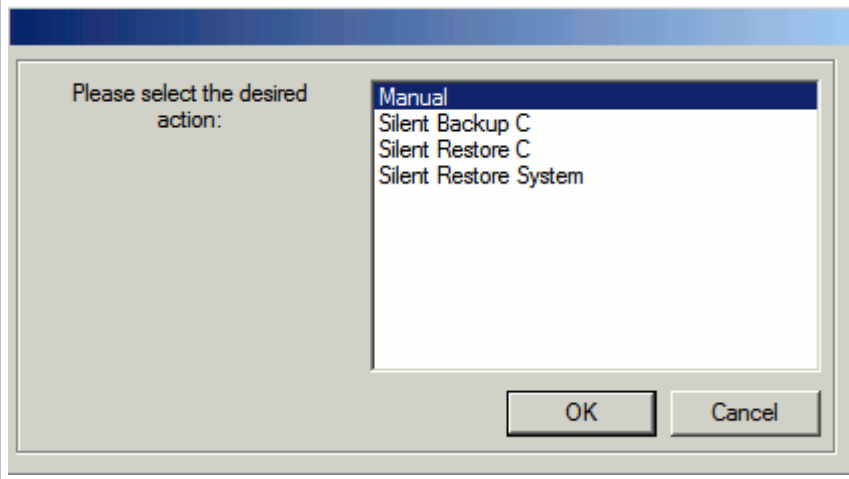
ACTION	MRAWP/MRWP
■ Wait until the COBR tool window is displayed. Fig. 18: COBR main Backup/Restore popup 	☐
■ Connect the USB medium.	☐
■ In the Backup/Restore pop-up select Manual and press OK Fig. 19: Manual Backup/Restore popup 	☐

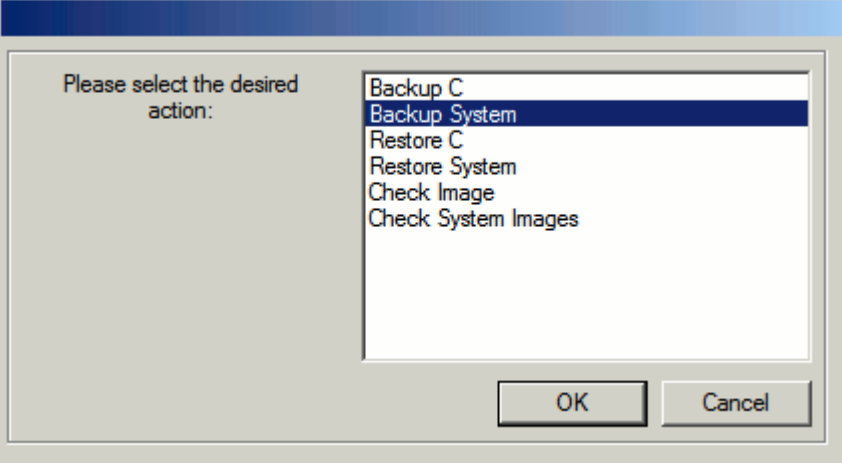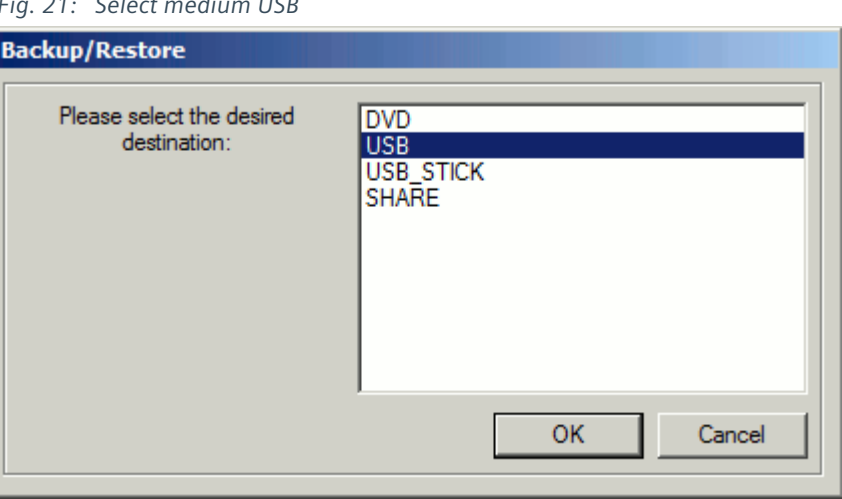ACTION	MRAWP/MRWP
In the next Backup/Restore pop-up select Backup System and press OK. Fig. 20: Backup/Restore popup - Backup System 	☐
In the next Backup/Restore popup select USB as the medium for the backup image and press OK. Fig. 21: Select medium USB 	☐

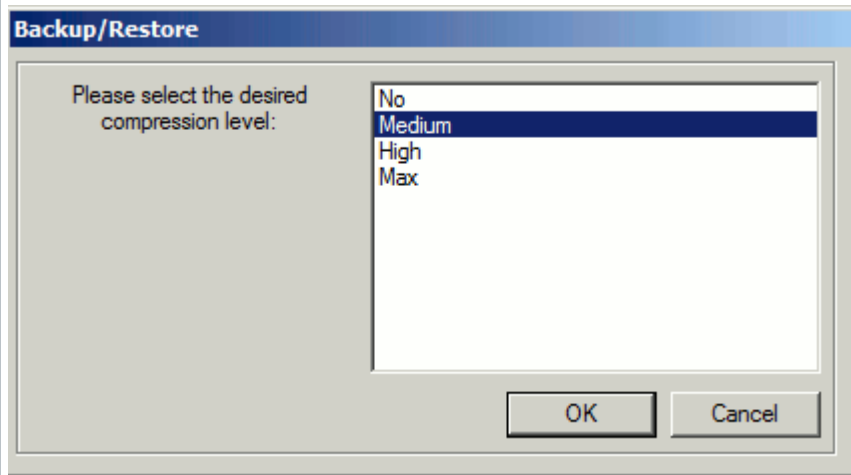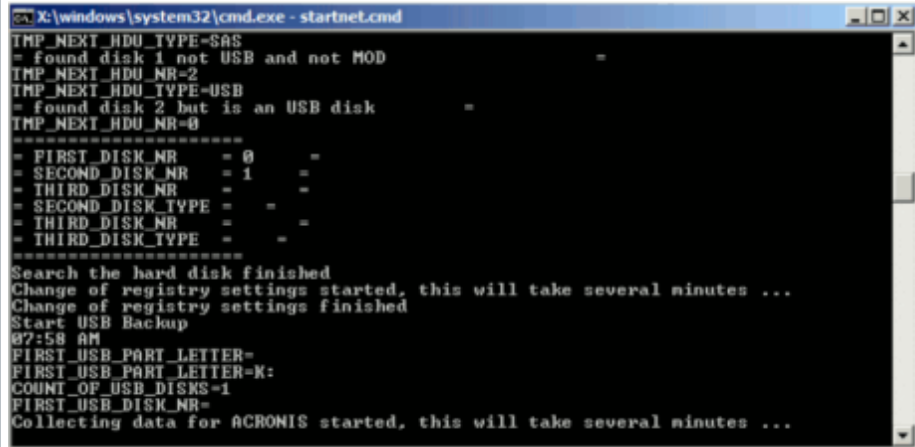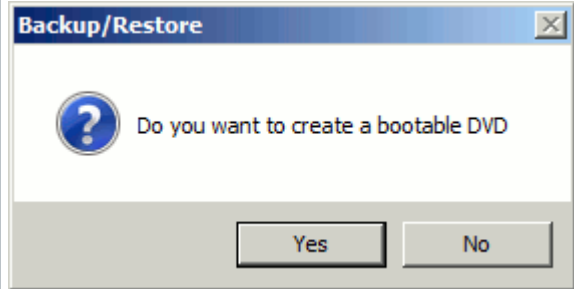ACTION	MRAWP/MRWP
In the next Backup/Restore pop-up select the compression level medium of the data and press OK. Fig. 22: Backup/Restore popup - compression level 	☐
Now the system back up will run. The following is an example of the screen display. Fig. 23: COBR Acronis image backup example output 	☐
When the backups of all disks are written answer the popup "Do you want to create a bootable DVD" with "No". Fig. 24: Create boot DVD pop-up 	☐

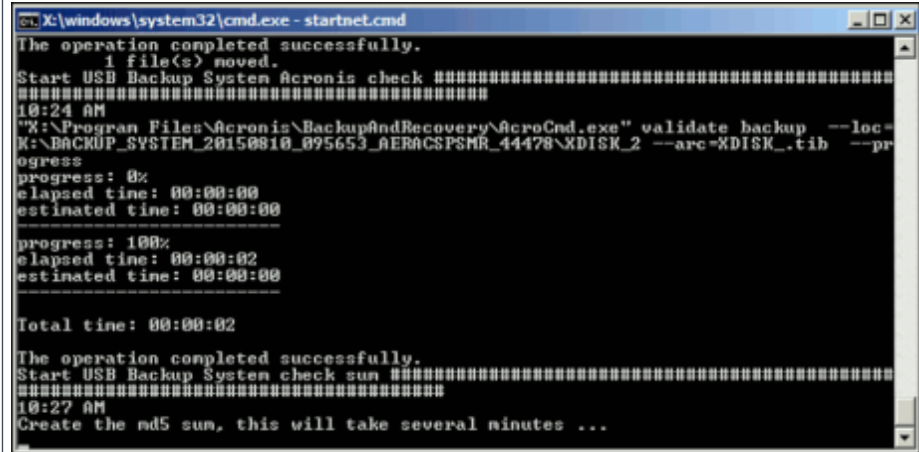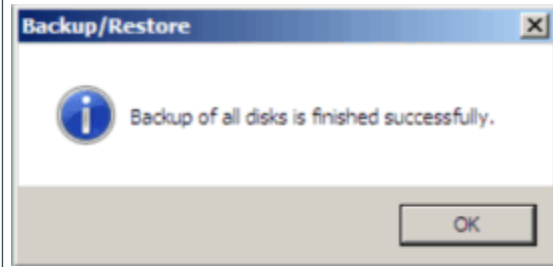ACTION	MRAWP/MRWP
Now a checksum (MD5) is created for each image file. This will take some minutes. Fig. 25: MD5 checksum creation 	☐
Confirm the success popup by pressing OK . Fig. 26: Backup finished successfully 	☐

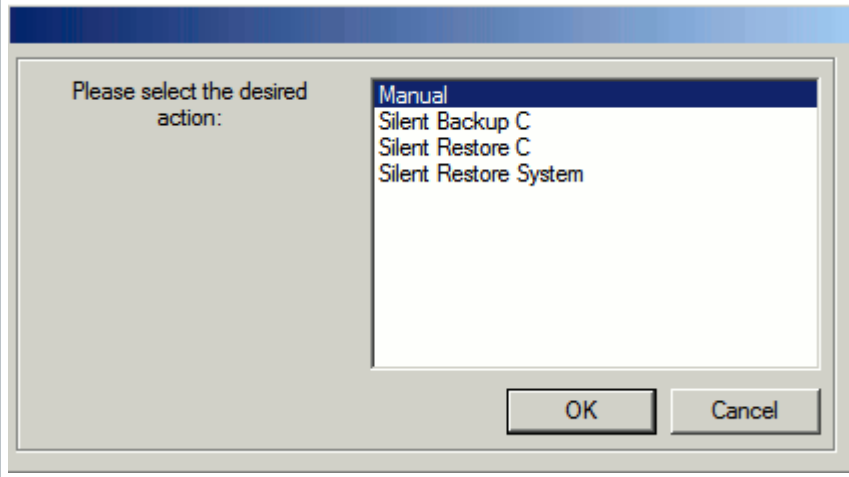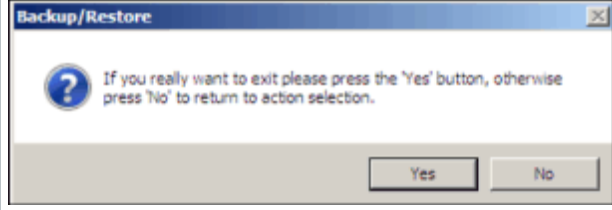ACTION	MRAWP/MRWP
■ The main Backup/Restore popup will appear. Fig. 27: COBR main Backup/Restore popup 	☐
■ Press Cancel to exit the Backup/Restore procedure. Fig. 28: Main Backup/Restore popup 	☐
■ Confirm the following popup with Yes. Fig. 29: Confirmation popup 	☐

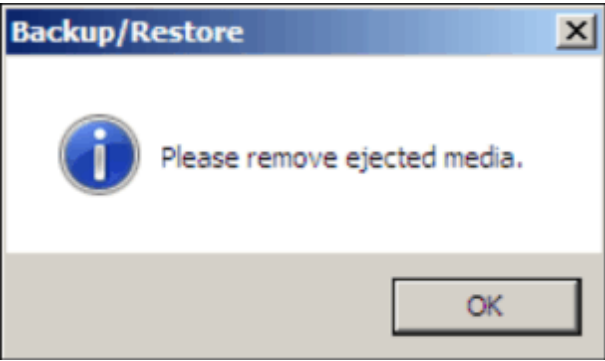ACTION	MRAWP/MRWP
The CD/DVD-ROM drive opens and the following popup is displayed. <i>Fig. 30: Backup/Restore remove ejected medium popup</i> 	☐
Remove the ejected disk from the drive.	☐
Close the drive and click OK to confirm.	☐
The system is rebooted.	☐
Disconnect the USB medium.	☐
After the reboot the windows login ist dispalyed. Please wait until a further reboot is performed (this takes about 30 seconds).	☐
The windows login is now displayed for the second time after the reboot. Now login as user meduser and wait until the main UI is displayed.	☐

3.1.3.3 Manual Backup of System Images to DVD


This procedure describes how to back up the complete software system for MR AWP (MRC) / MR WP (MRSC) (all hard disks and partitions) to DVD medias.




Prerequisites



- DVD's (about 10)
- Autologin is disabled
- Masterbackup is created and stored to partition ASIS drive letter **F:**.
- The step (→ File cleanup for image creation / Page 128) is performed


Procedure



ACTION	MRAWP/MRWP
Reboot and log in as user administrator	☐
Insert the first DVD from the last created backup image into the CD/DVD-ROM drive R:	☐


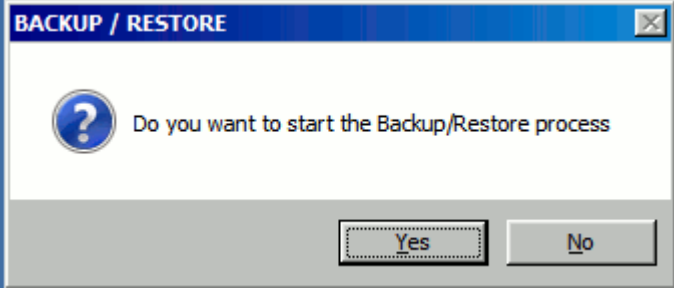
ACTION	MRAWP/MRWP
■ Invoke a command shell (Start ⇒ All Programs ⇒ Accessories ⇒ Command Prompt)	☐
■ In the command shell, type in: <code>C:\MedCom\bin\StartDiskImageGeneration.exe</code> ↵ The command above will prepare the system for COBR image creation.	☐
■ The COBR tool DVD is requested. Insert the first DVD from the last created backup image into the CD/DVD-ROM drive.	☐
■ Press any key e.g. <SPACE> to start the COBR image backup.	☐
■ Confirm the displayed popup with Yes. <i>Fig. 31: Backup/Restore popup</i> 	☐
■ The running script performs a reboot at the end. This may take some time. Please wait!.	☐
■ A reboot is performed and the COBR tool boots.	☐

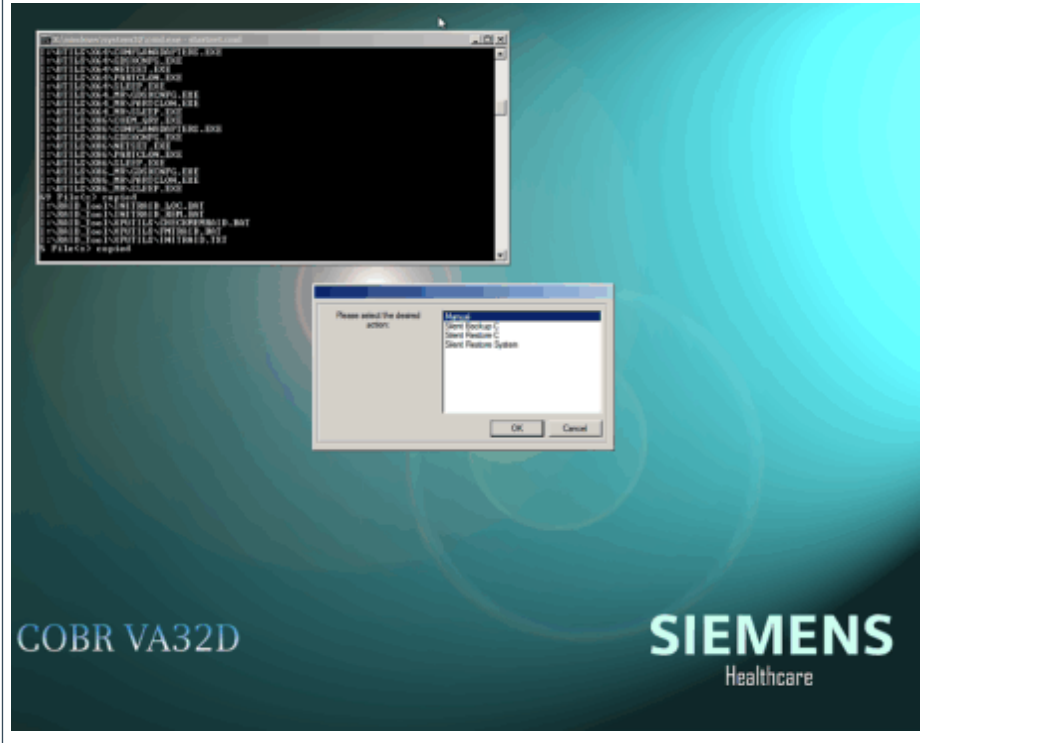
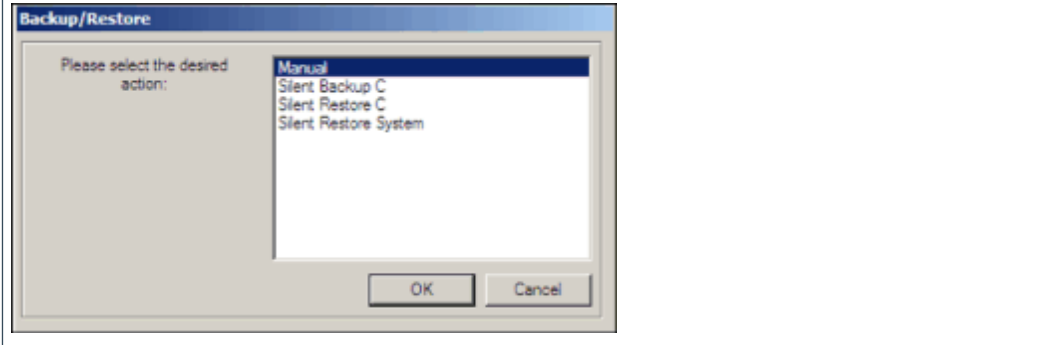
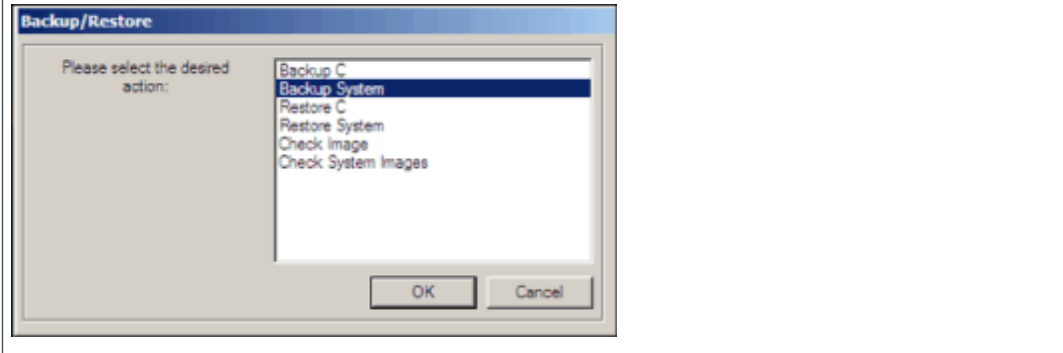
ACTION	MRAWP/MRWP
Wait until the COBR tool window is displayed. Fig. 32: COBR main Backup/Restore popup 	☐
In the Backup/Restore pop-up select Manual and press OK Fig. 33: Manual Backup/Restore popup 	☐
In the next Backup/Restore pop-up select Backup System and press OK. Fig. 34: Backup/Restore popup - Backup System 	☐

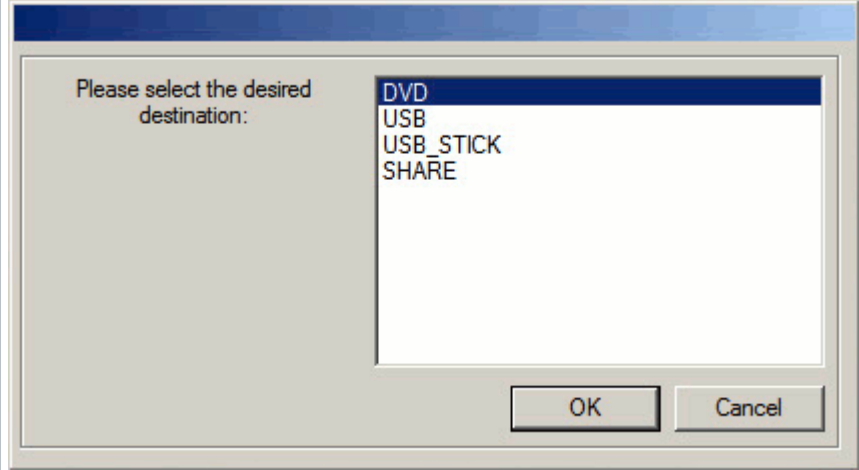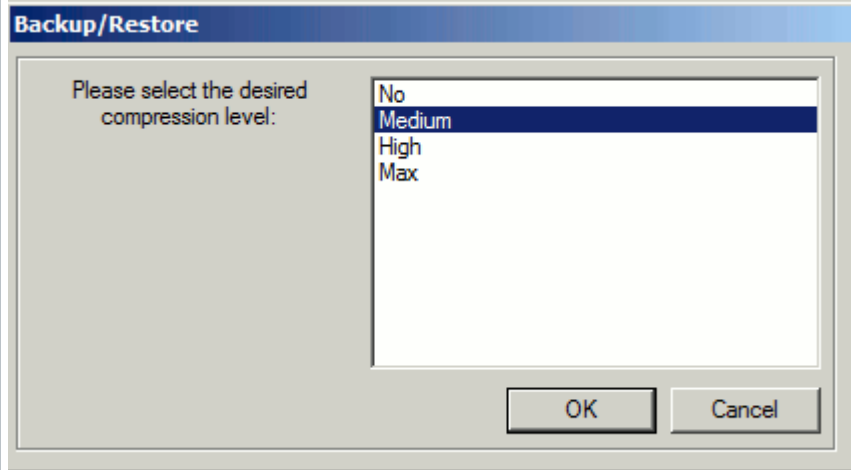
ACTION	MRAWP/MRWP
■ In the next Backup/Restore popup select DVD as the medium for the backup image and press OK. Fig. 35: Select medium DVD 	☐
■ In the next Backup/Restore pop-up select the compression level medium of the data and press OK. Fig. 36: Backup/Restore popup - compression level 	☐

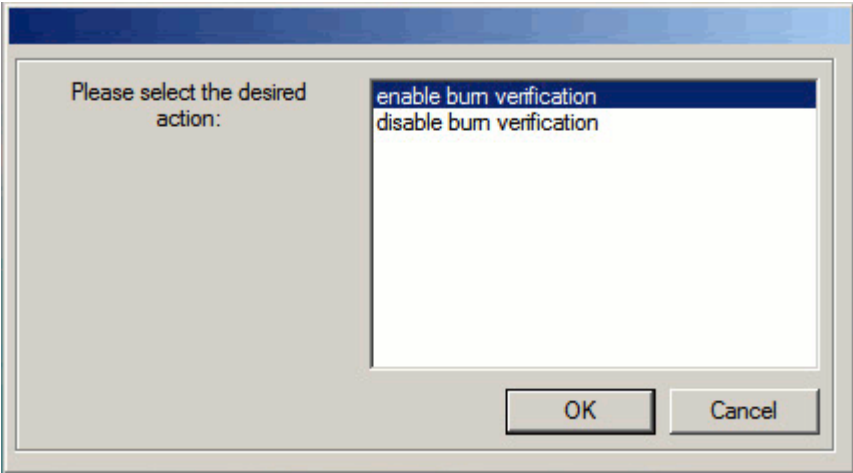ACTION	MRAWP/MRWP
■ In the next Backup/Restore pop-up select enable burn verification and press OK. Fig. 37: COBR backup enable burn verification 	☐
■ Now the system back up will run. The following is an example of the screen display. Fig. 38: COBR Acronis image backup example output 	☐

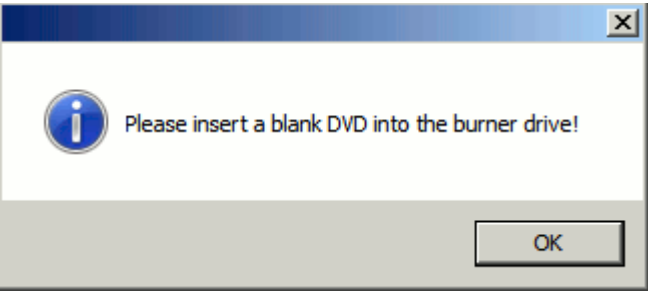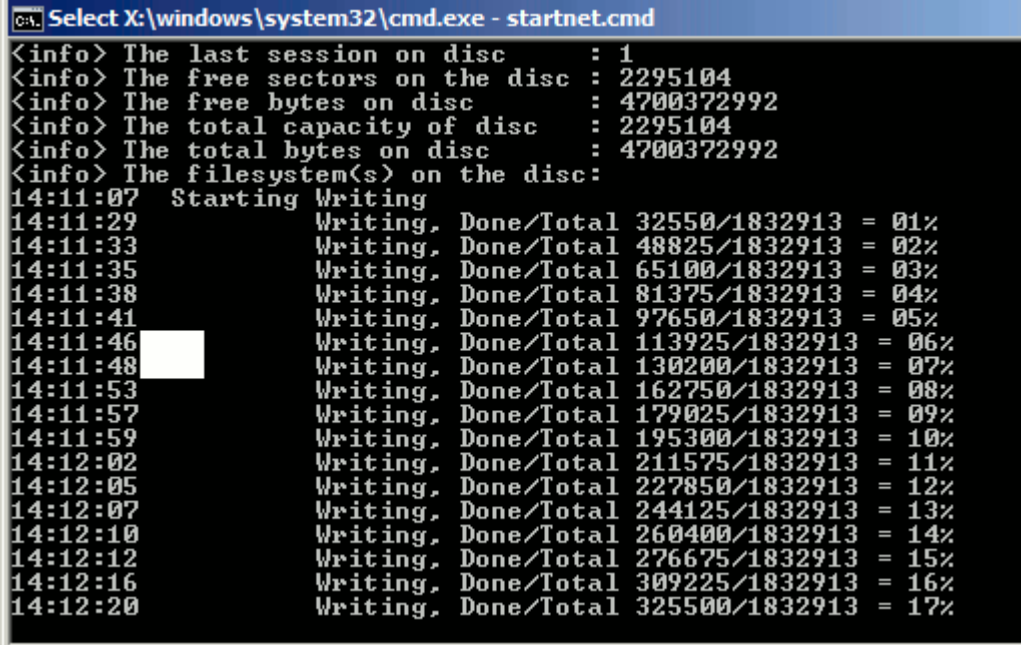ACTION	MRAWP/MRWP
Wait until the following request for an empty DVD is displayed: <i>Fig. 39: COBR screen output - request DVD</i> 	☐
- Label this first DVD as follows: - System serial number + Date + Softwareversion + System Image + DVD #1 of n. n is the number of DVDs used in the image DVD set. Write down this number at the end when all DVDs have been created	☐
The CD/DVD-R drive opens. Insert this labeled empty DVD into this drive and close the drive.	☐
Confirm the displayed popup with OK.	☐
DVD generation now starts. The following shows an example of it. <i>Fig. 40: DVD generation - writing data to DVD</i> 	☐

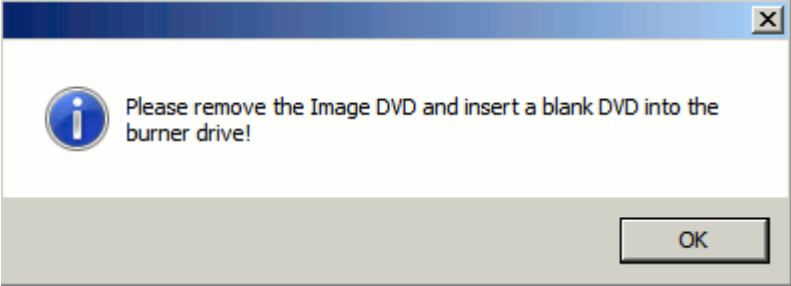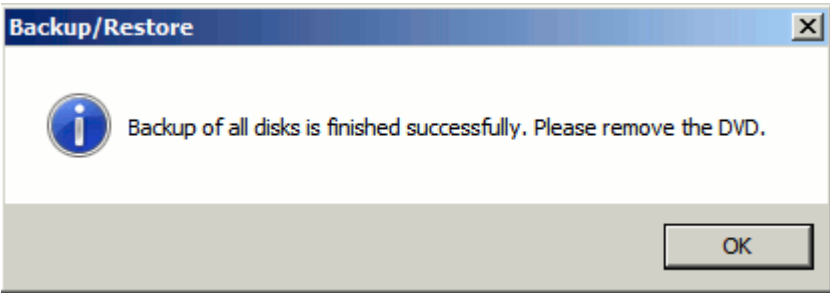ACTION	MRAWP/MRWP
An additional DVD is requested automatically if needed. The following popup shows the request. <i>Fig. 41: DVD generation - request popup</i> 	☐
- Label the additional DVDs as follows: - Second DVD: System serial number + Date + Softwareversion + System Image + DVD #2 + n - Third DVD if needed: System serial number + Date + Softwareversion + System Image + DVD #3 + n - Any additional DVDs if needed: System serial number + Date + Softwareversion + System Image + DVD #x + n x is the number starting with 4 and incremented by 1 with each additional DVD request. n is the number of DVDs used in the image DVD set. Write down this number at the end when all DVDs have been created	
Remove the ejected DVD, insert an empty DVD into the drive, and close the drive. Confirm this popup with OK. Continue with this DVD as with the first one. Repeat this procedure until all data are written to the DVDs.	☐
Once all data are successfully written to DVD the following popup is displayed. <i>Fig. 42: Backup/Restore to DVD successfully finished popup</i> 	☐
Remove the ejected medium, close the drive, and confirm this popup with OK.	☐

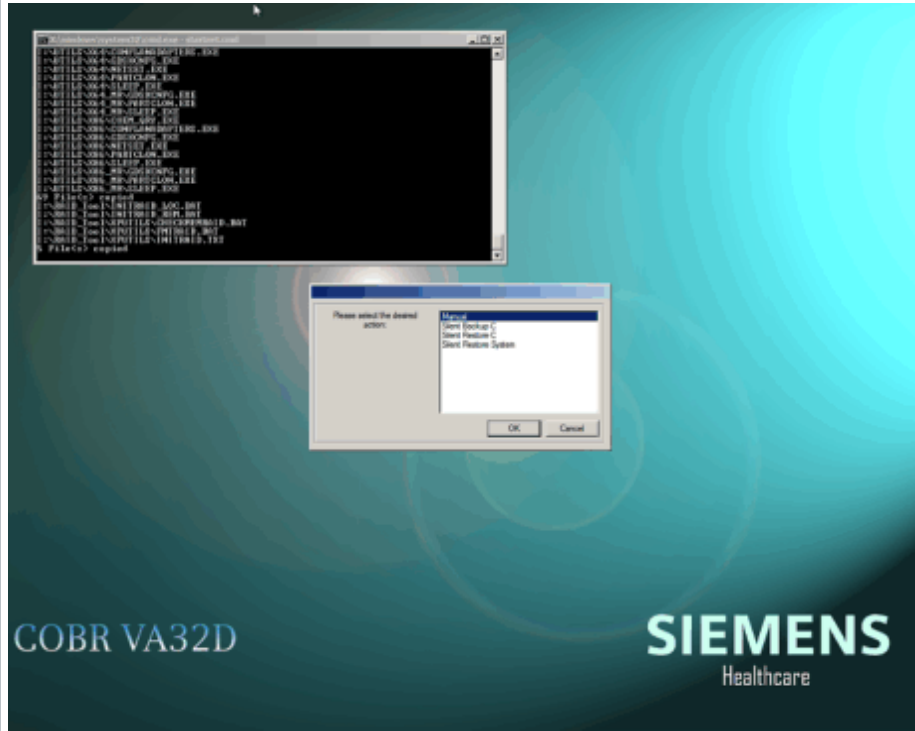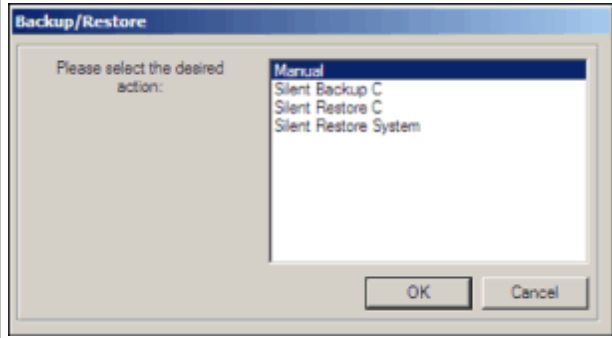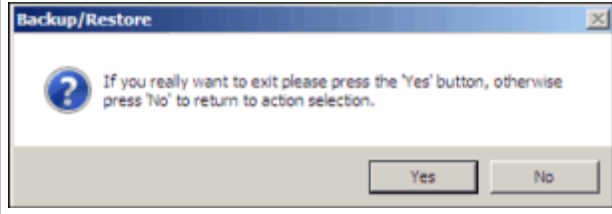ACTION	MRAWP/MRWP
■ The main Backup/Restore popup will appear. Fig. 43: COBR main Backup/Restore popup The screenshot shows the COBR VA32D Siemens Healthcare interface. In the background, there is a command window with a list of commands and their outputs. In the foreground, a 'Backup/Restore' popup is displayed, asking the user to select an action from a list: Manual, Silent Backup C, Silent Restore C, and Silent Restore System. The popup has 'OK' and 'Cancel' buttons.	☐
■ Press Cancel to exit the Backup/Restore procedure. Fig. 44: Main Backup/Restore popup The screenshot shows the 'Backup/Restore' popup with the title bar. The text 'Please select the desired action:' is followed by a list box containing 'Manual', 'Silent Backup C', 'Silent Restore C', and 'Silent Restore System'. At the bottom are 'OK' and 'Cancel' buttons.	☐
■ Confirm the following popup with Yes. Fig. 45: Confirmation popup The screenshot shows a 'Backup/Restore' confirmation popup. It contains a question mark icon and the text: 'If you really want to exit please press the 'Yes' button, otherwise press 'No' to return to action selection.' At the bottom are 'Yes' and 'No' buttons.	☐

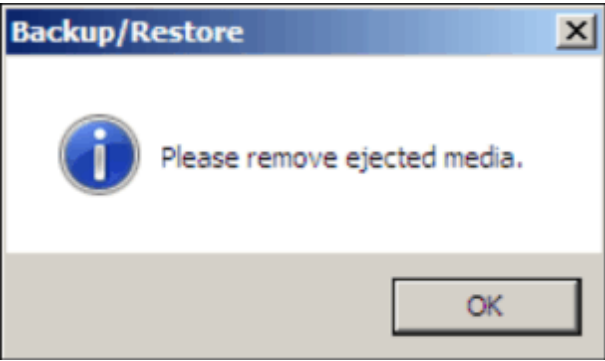ACTION	MRAWP/MRWP
The CD/DVD-ROM drive opens and the following popup is displayed. <i>Fig. 46: Backup/Restore remove ejected medium popup</i> 	☐
Remove the ejected disk from the drive.	☐
Close the drive and click OK to confirm.	☐
The system is rebooted.	☐
After the reboot the windows login ist dispalyed. Please wait until a further reboot is performed (this takes about 30 seconds).	☐
The windows login is now displayed for the second time after the reboot. Now login as user meduser and wait until the main UI is displayed.	☐

3.1.4 Restore of the COBR Image

3.1.4.1 Manual Restore of Partition C Image from USB-Disk

This procedure describes how to restore a complete software image backup or MR AWP (MRC) / MR WP (MRSC) (all hard disks and partitions) from a USB-disk.



Prerequisites




- USB-disk with actual COBR image backup
- No login is possible (softwer is erased with the Blancco tool or spare part host)
- MR system is switched off at the alarmbox



Procedure no login as sysadmin is possible




- COBR image restore steps



ACTION	MRAWP/MRWP
Remove the power line connector from the MRIR in the Electronic Power Cabinet.	☐
Switch on the MR system at the alarmbox .	☐

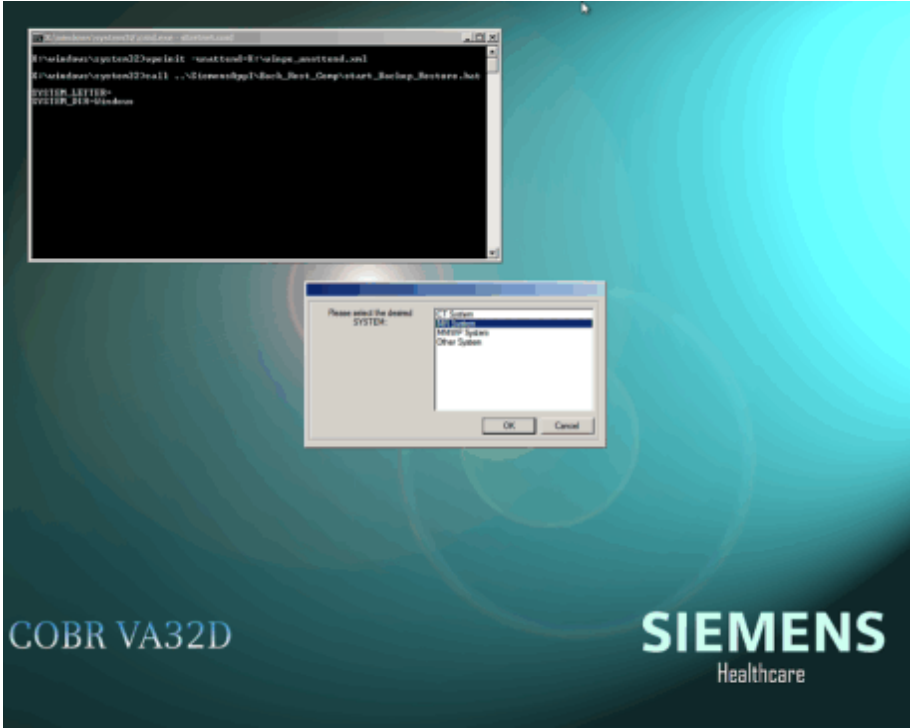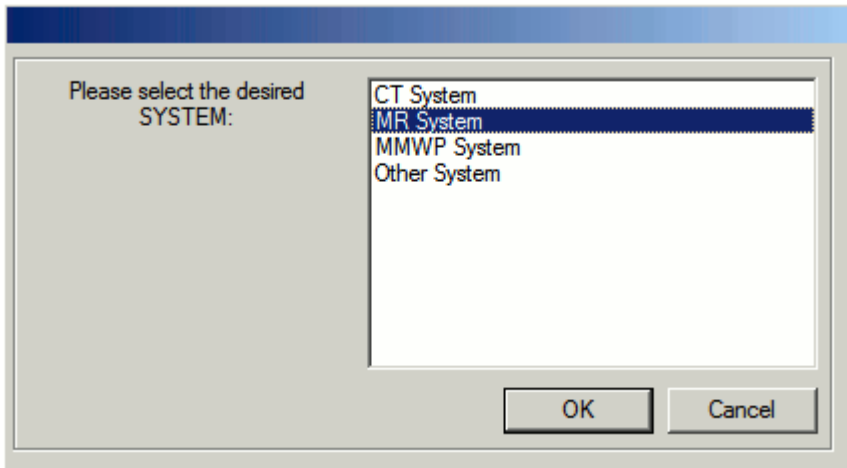ACTION	MRAWP/MRWP
Insert the first DVD (includes the COBR tool) from the last created backup image (DVD version) or the COBR tool DVD into the CD/DVD-ROM drive.	☐
The COBR tool boots.	☐
Wait until the COBR tool window is displayed. Fig. 47: COBR main Backup/Restore popup 	☐
In the displayed popup the SYSTEM modality is requested. Select MR System and press OK Fig. 48: Modality request popup 	☐
Connect the USB medium which contains the COBR image.	☐

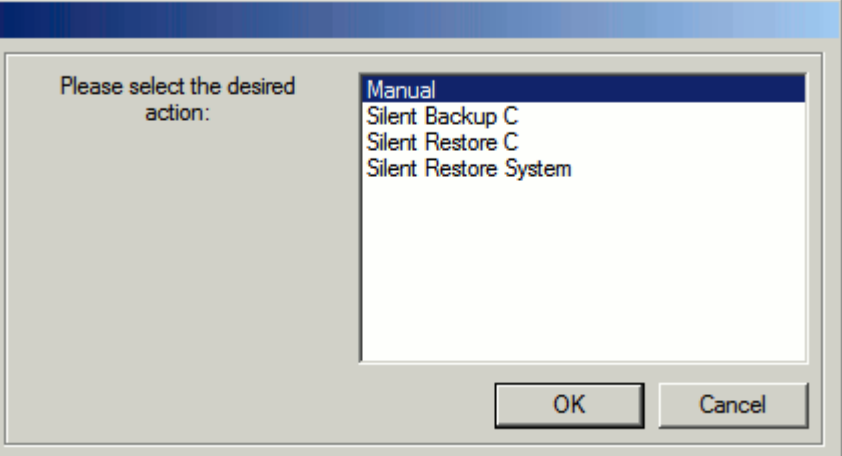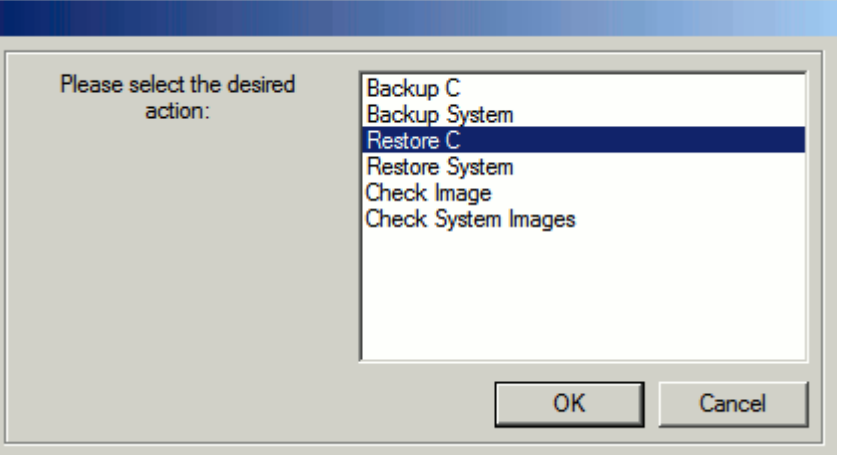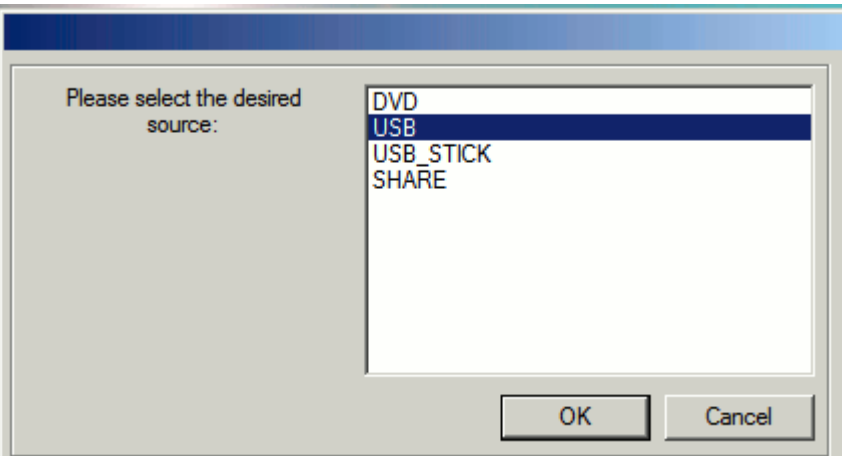ACTION	MRAWP/MRWP
In the next displayed pop-up select Manual and press OK <i>Fig. 49: Manual popup</i> 	☐
In the next pop-up select Restore System and press OK. <i>Fig. 50: Popup - Restore C</i> 	☐
In the next Backup/Restore pop-up select USB and press OK. <i>Fig. 51: Popup - select USB as source</i> 	☐

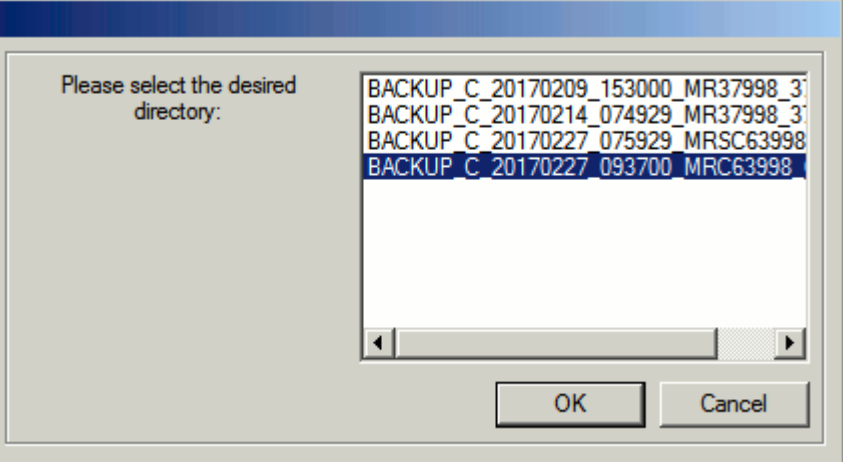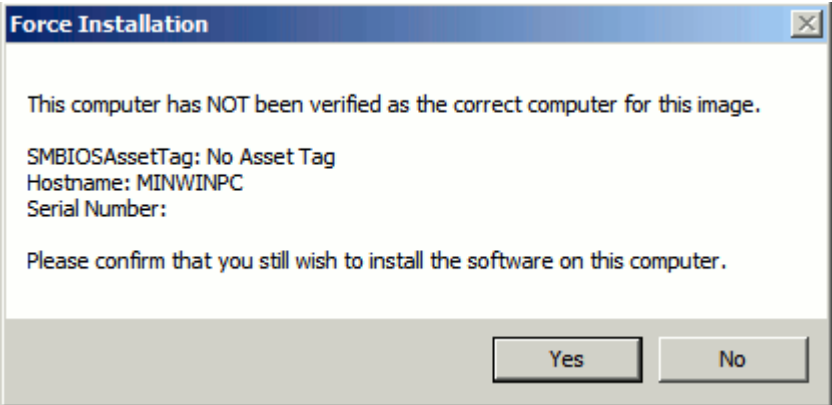ACTION	MRAWP/MRWP
In the next popup select the COBR directory which you want to restore (this directory contains the desired images) and press OK. Fig. 52: Example for the image restore selection 	☐
Confirm the Force Installation popup (example is shown below) with Yes. Fig. 53: Example of the force Installation popup  This Force Installation popup is displayed twice. Confirm it each time with Yes.	☐

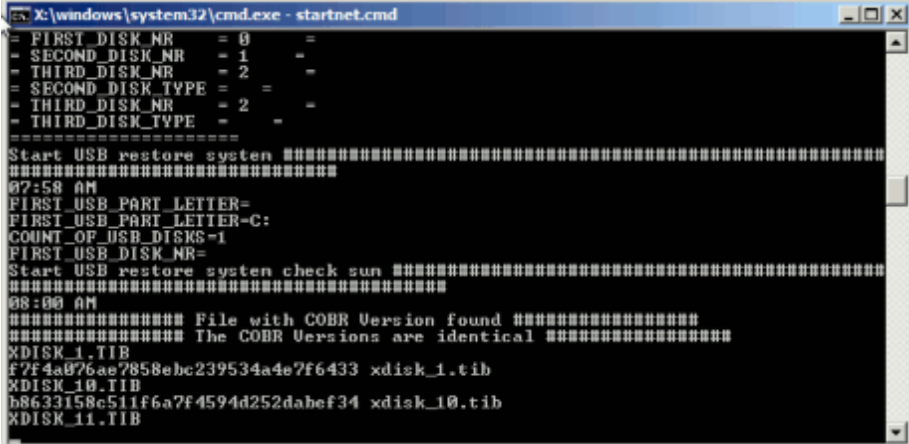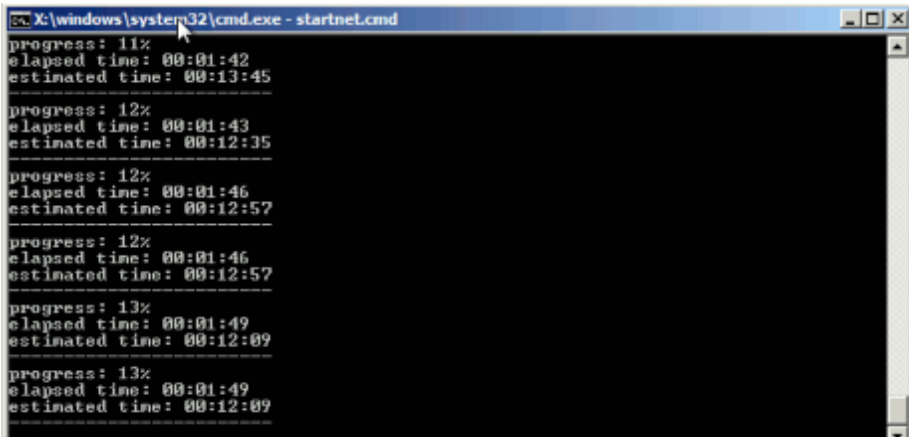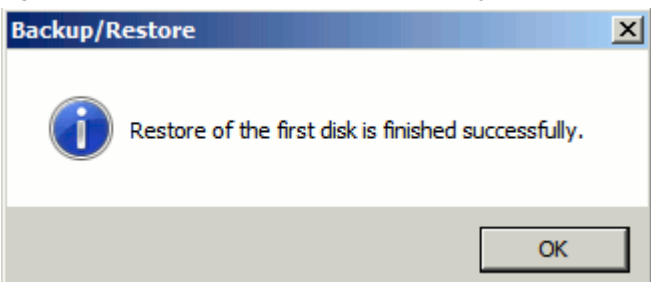ACTION	MRAWP/MRWP
Now the system restore will run. First the md5 verification is performedThe following is an example of the screen display. Fig. 54: COBR md5 check example output  Fig. 55: COBR image restore example output 	☐
Confirm the success popup by pressing OK . Fig. 56: Restore first disk finished successfully 	☐

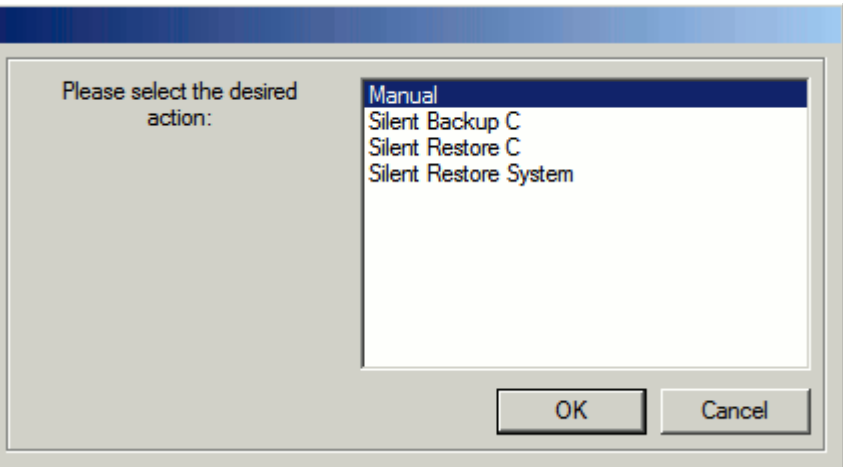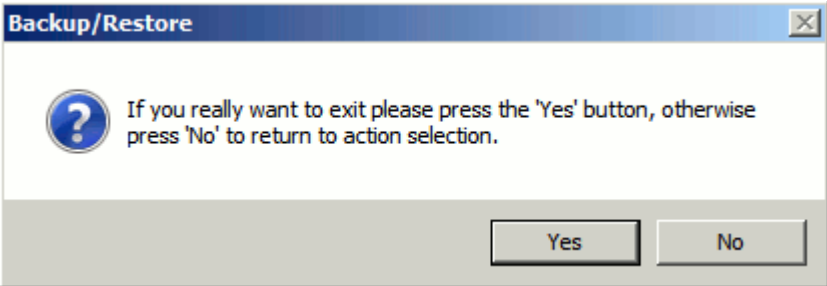ACTION	MRAWP/MRWP
The main popup will appear. <i>Fig. 57: COBR main popup</i> 	☐
Press Cancel to exit the Backup/Restore procedure.	☐
Confirm the following popup with Yes. <i>Fig. 58: Confirmation popup</i> 	☐
The CD/DVD-ROM drive opens and the following popup is displayed. <i>Fig. 59: Backup/Restore remove ejected medium popup</i> 	☐
Remove the ejected disk from the drive.	☐
Close the drive and click OK to confirm.	☐
The system is rebooted.	☐
Disconnect the USB medium.	☐

- COBR post restore steps on MRAWP (MRC)


ACTION	MRAWP
Note: The following behaviour is new in the VB19B COBR image restore.	
■ After the reboot, the autologin fails with an error message. Please wait until a further reboot is done (taks about 30 seconds). The next autologin is successful done. Note: Please don't klick or hit any key ! Just wait.	☐
■ After the autologin, additional configuration steps are performed. Wait until they are finished, and the system is rebooted. This will take about 10 minutes. Note: Don't accept any displayed popup or press a button on any displayed popup.	☐

ACTION

MRAWP


During these steps, do not close any displayed command windows. Also do not confirm any displayed popup windows. The following windows will be displayed:

Fig. 60: Displayed popup

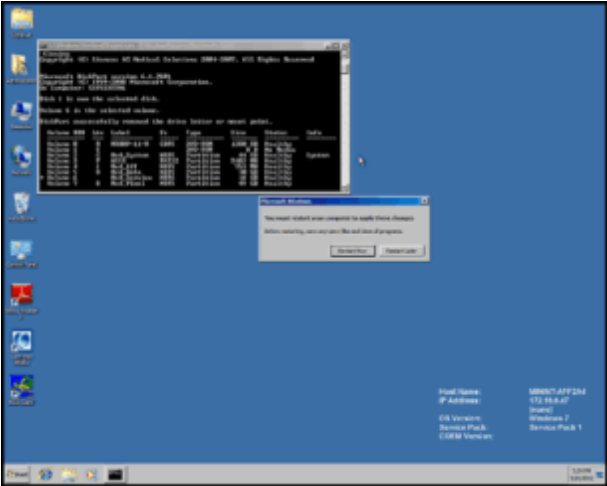


Note: Don't accept any displayed popup or press a button on any displayed popup. If the following popup is displayed, don't touch or accept it. Keep it displayed on the screen. It will be closed automatically.

Fig. 61: Windows popup

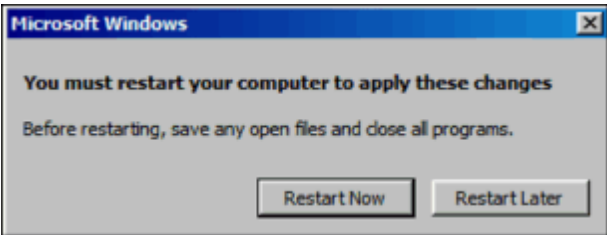

Note: The following popup can be ignored too, if it is displayed.

Fig. 62: Network location popup

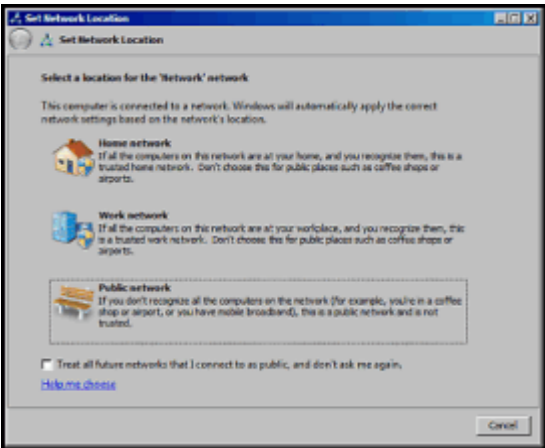


ACTION

MRAWP

The CSSI Post popup is displayed (see) (→ Fig. 63 Page 45)



Fig. 63: CSSI Postaction popup

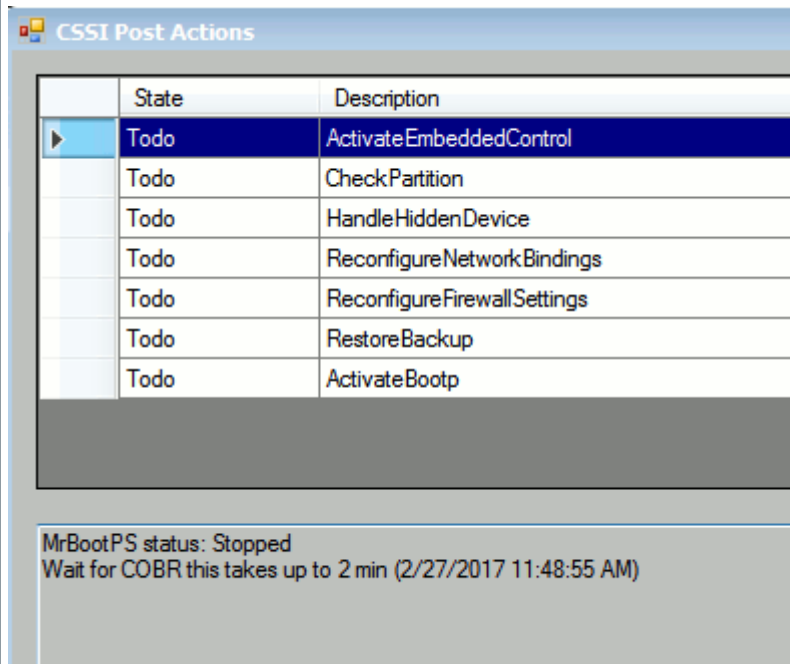


The following CSSI Post Action steps are performed:

- ActivateEmbeddedControl (duration under 1 min.)
- CheckPartition (duration under 1 min.)

This check will fail, because only the partition C is available.

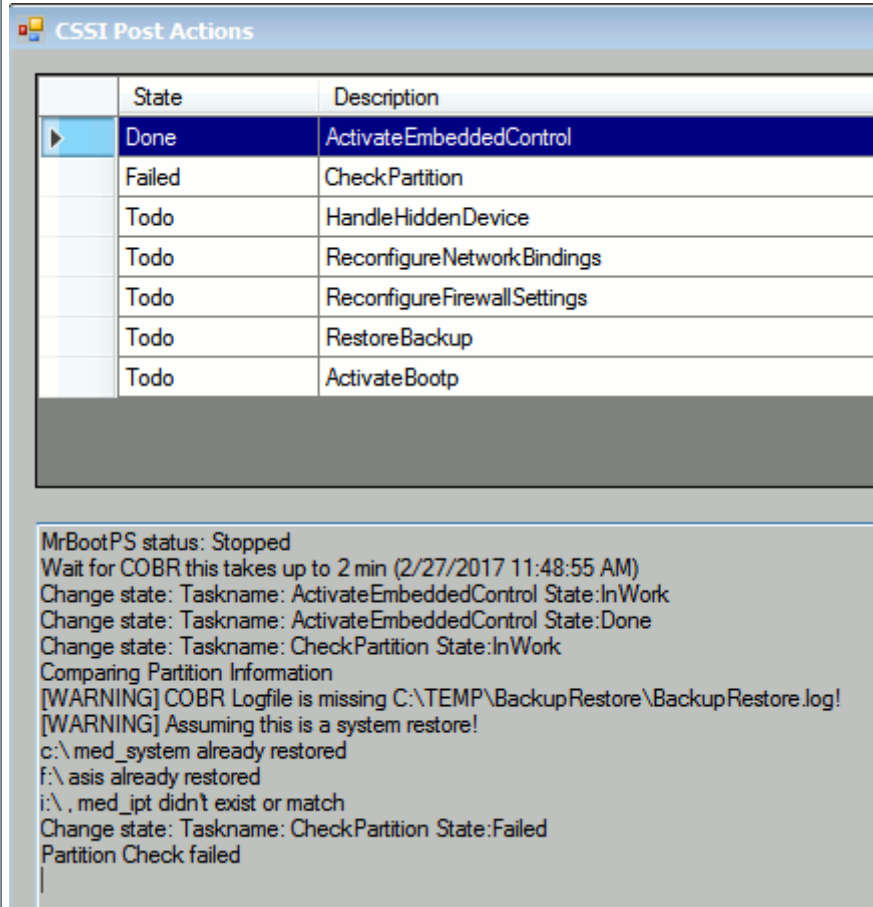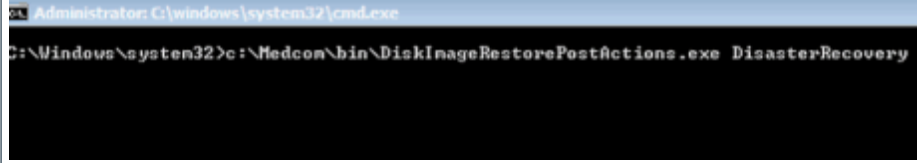ACTION	MRAWP
Please wait until you get the following message in the CSSI Post Action Popup: Fig. 64: CSSI Post Action Partition Check fails 	☐
Close the CSSI Post Actions by pressing the X.	☐
Select one of the displayed cmd tool windows.	☐
Type in the following comand: C:\MedCom\bin\DiskImageRestorePostActions.exe DisasterRecovery Fig. 65: Input command 	☐

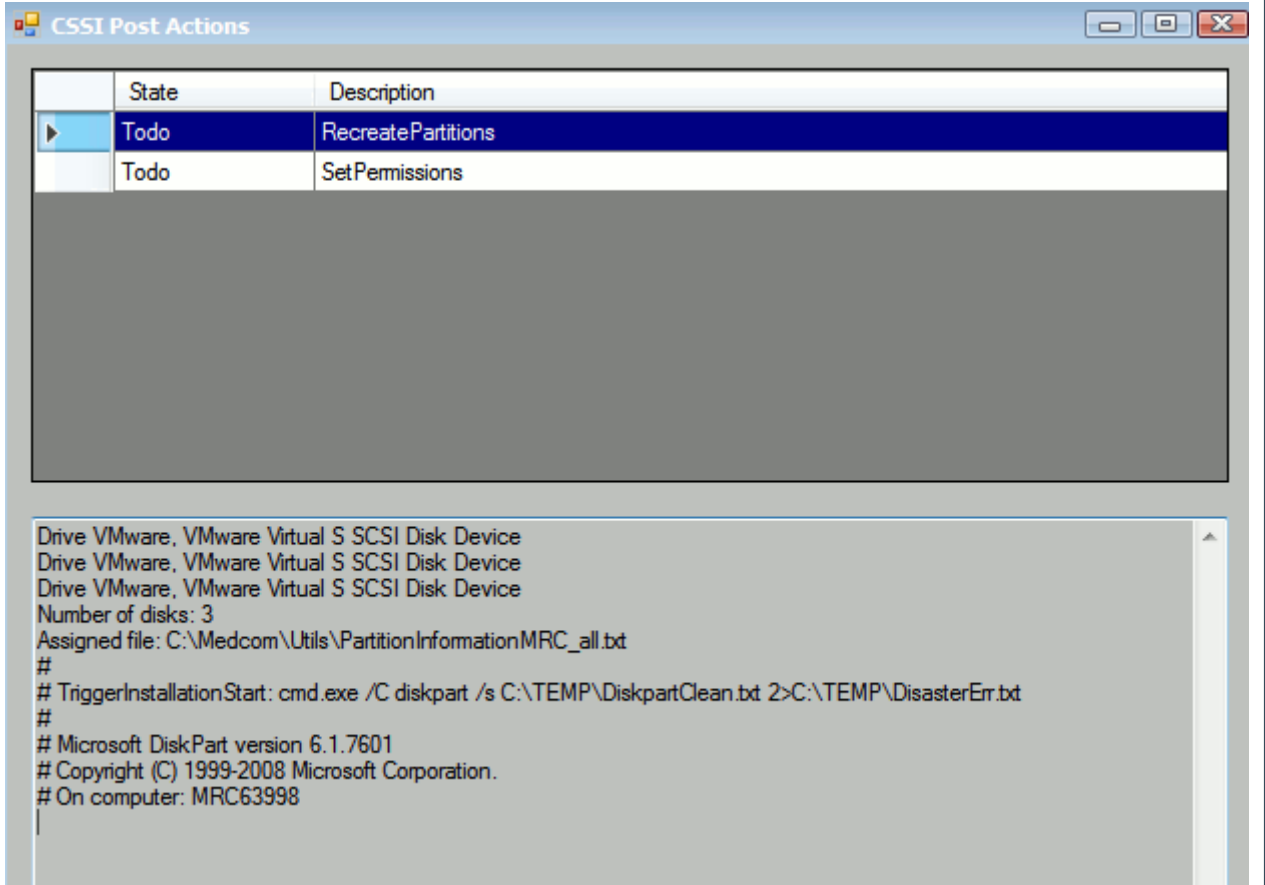ACTION		MRAWP
After Disaster Recovery tool is launched, the GUI with following tasks is displayed on the screen:		
Fig. 66: DesasterRecovery Popup		
		
■ The missing partitions are created and formatted and the databases are recreated automatically. No user interaction is required. Two reboots are done.		☐
■ A autologin as meduser is performed after the second reboot. Wait until the main UI is up.		☐
■ Connect the power line connector from the MRIR in the Electronic Power Cabinet		☐

ACTION	MRAWP
Wait until the measurement is ready.	☐
Note: In case of an error the CSSI Post Action stops at the step which shows failed as status. The reboot is than not done automatically, please perform the following workaround in this case. Run the CSSI Post Actions a second time. - Close the CSSI Post Actions by pressing the X. - Start a cmd tool window (Start ⇒ All Programs⇒ Command Prompt). - In the cmd tool window change the directory to C:\MedCom\bin - Execute again the CSSI Post actions by typing in the command: <code>DiskImageRestorePostActions.exe ↵</code> Now the CSSI Post actions starts again. This takes about 20 minutes. After this time a reboot is done automatically. If the CSSI Post actions tool fails again, please contact the HSC MR in Erlangen.	

- COBR post restore steps on MRWP (MRSC)


ACTION	MRWP
Note: The following behaviour is new in the VB19B COBR image restore.	
After the reboot, the autologin fails with an error message. Please wait until a further reboot is done (taks about 1 minute) The next autologin is successful done. Note: Please don't klick or hit any key ! Just wait.	☐
After the autologin, additional configuration steps are performed. Wait until they are finished, and the system is rebooted. This will take about 5 minutes. Note: Don't close or confirm any displayed popups.	☐

ACTION

MRWP

During these steps, do not close any displayed command windows. Also do not confirm any displayed popup windows. The following windows will be displayed:


Fig. 67: Displayed popup

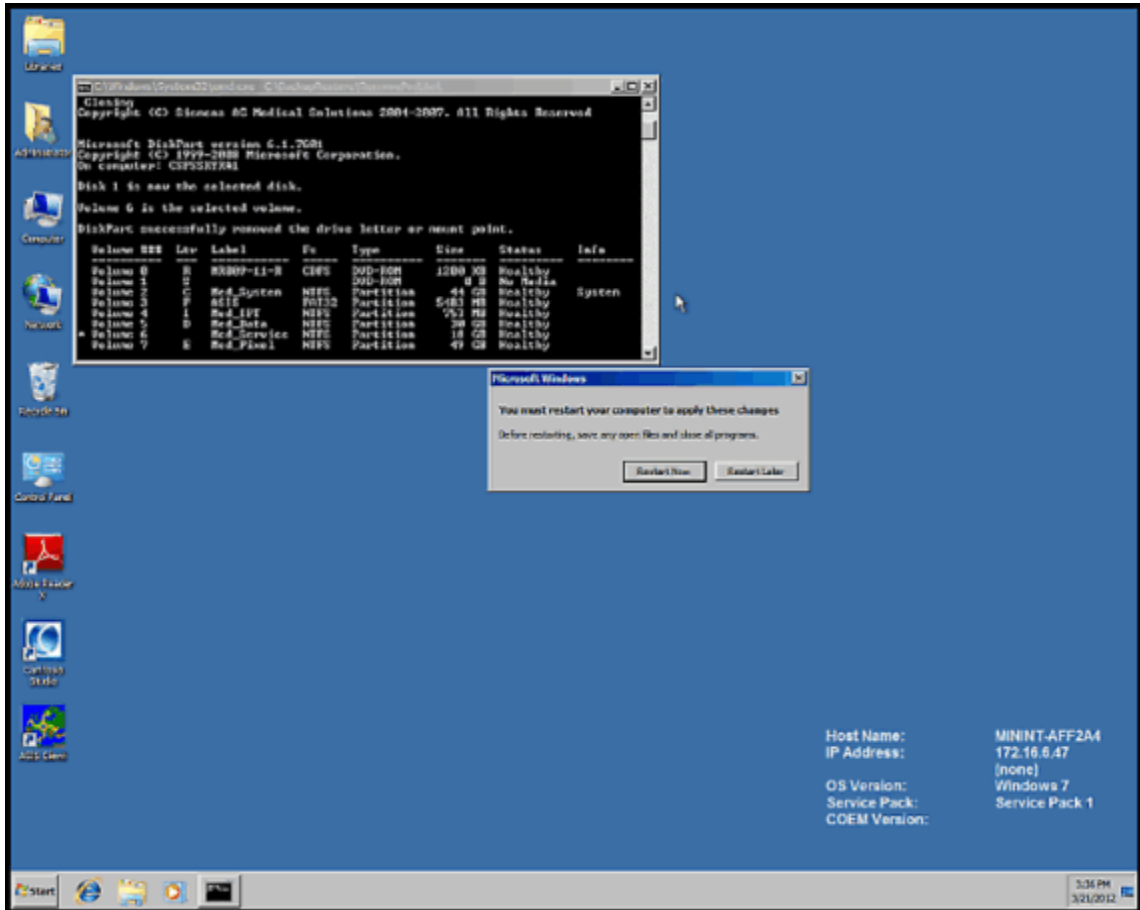

The CSSI Post popup is also displayed (see Fig. 60 “CSSI Post Actions popup”)

Note: Don’t accept any displayed popup or press a button on any displayed popup. If the following popup is displayed, don’t touch or accept it. Keep it displayed on the screen. It will be closed automatically.


Fig. 68: Windows popup

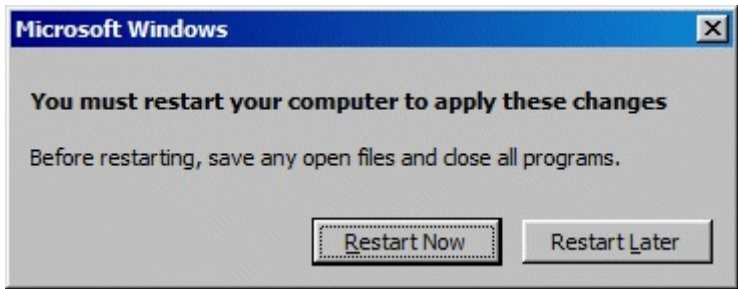

ACTION


MRWP

Note: The following popup can be ignored too, if it is displayed.
Fig. 69: Network location popup

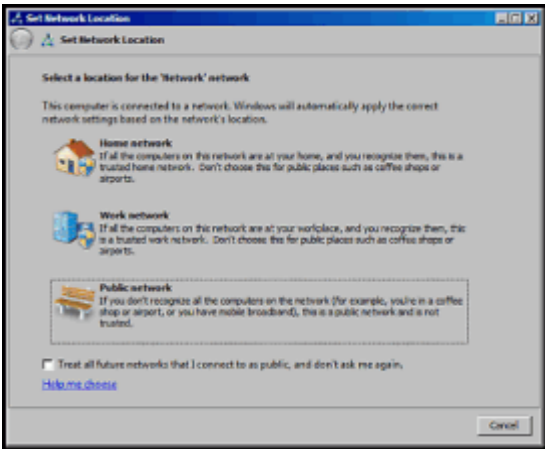

- The system is rebooted. ☐
- A autologin as user "syadmin" is performed and the following CSSI post popup is displayed. ☐



The CSSI Post popup is displayed (see)(→ Fig. 63 Page 45)
Fig. 70: CSSI Postaction popup

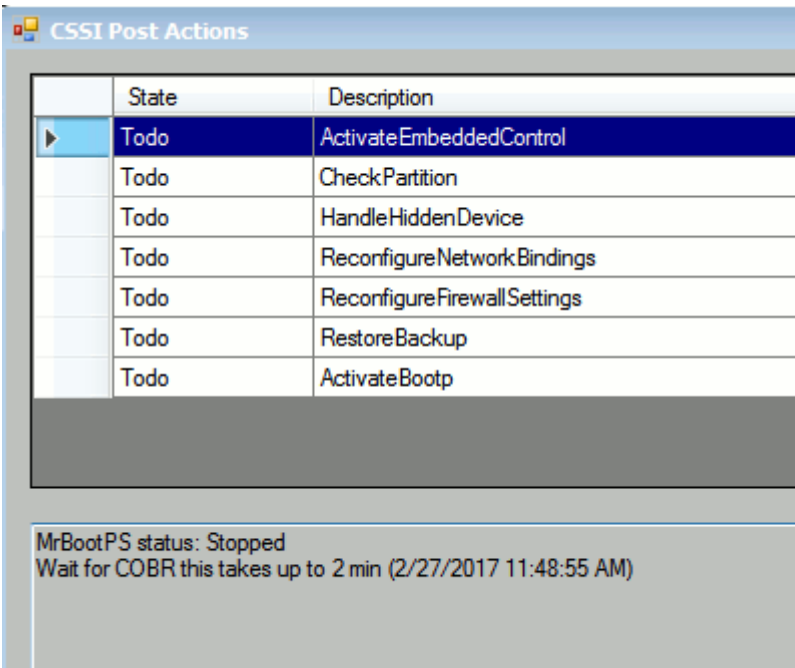

The following CSSI Post Action steps are performed:

- ActivateEmbeddedControl (duration under 1 min.)
- CheckPartition (duration under 1 min.)


This check will fail, because only the partition C is available.

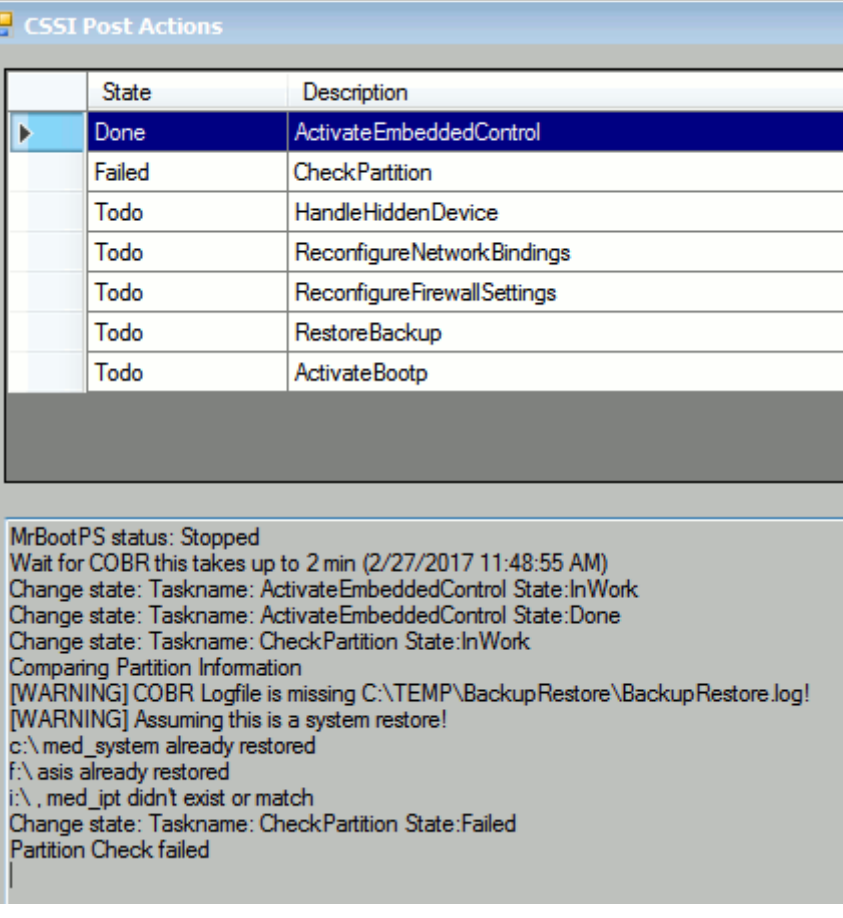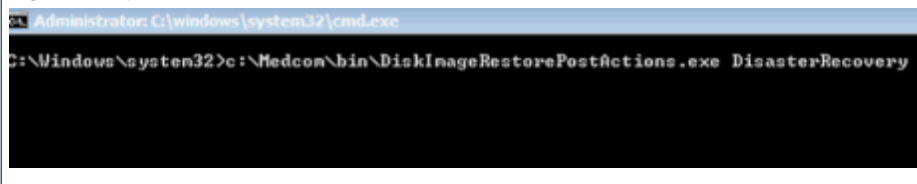ACTION	MRWP
Please wait until you get the following message in the CSSI Post Action Popup: <i>Fig. 71: CSSI Post Action Partition Check fails</i> 	☐
Close the CSSI Post Actions by pressing the X.	☐
Start a cmd tool windows from the Start menu.	☐
Type in the following comand: <code>C:\MedCom\bin\DiskImageRestorePostActions DisasterRecovery</code> <i>Fig. 72: Input command</i> 	☐

ACTION

MRWP

After Disaster Recovery tool is launched, the GUI with following tasks is displayed on the screen:
Fig. 73: DesasterRecovery Popup MRWP

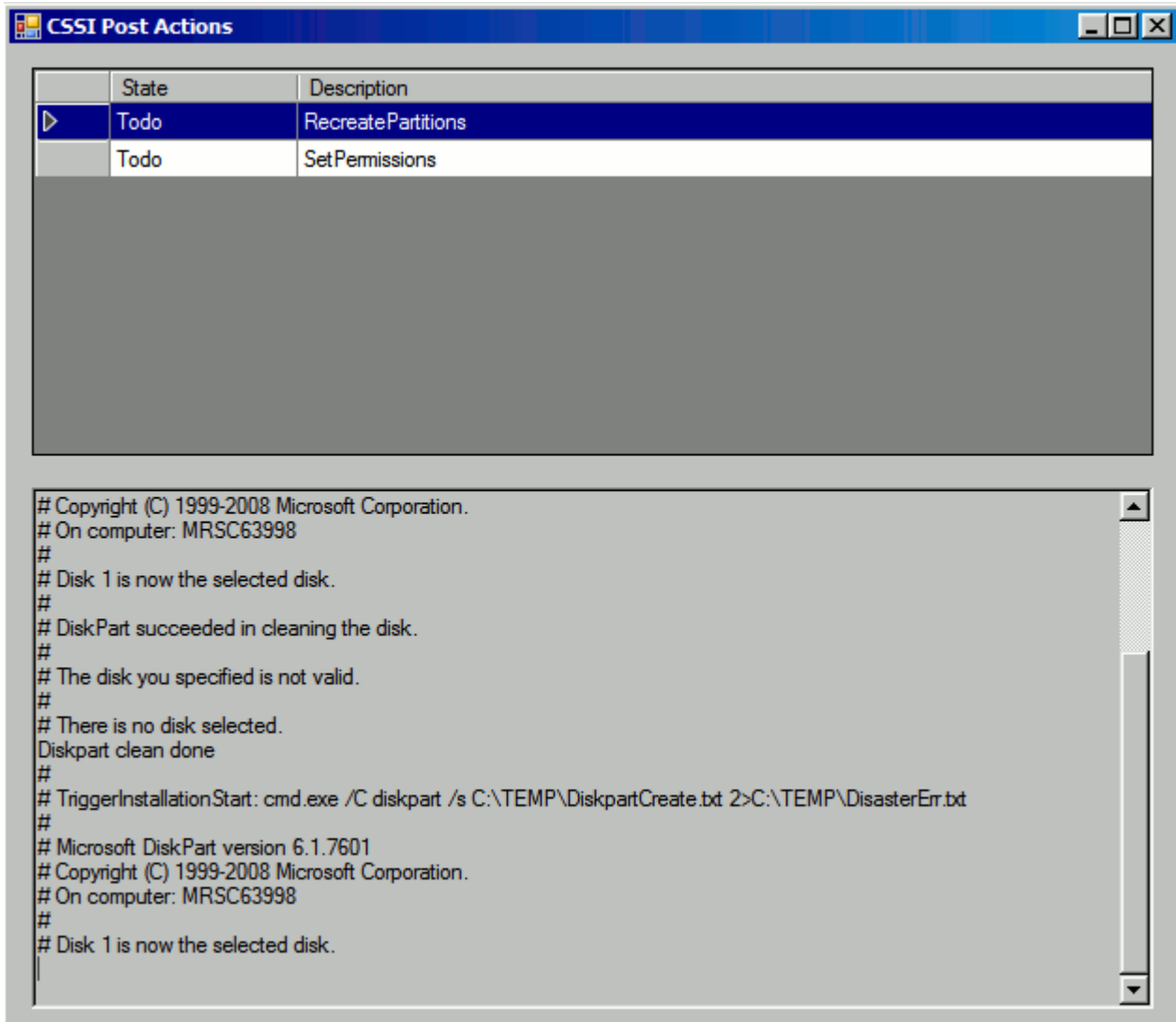

- The missing partitions are created and formatted and the databases are recreated automatically. No user interaction is required. Two reboots are done.
- ☐

ACTION	MRWP
■ A autologin as meduser is performed after the second reboot. Wait until the main UI is up.	☐
Note: In case of an error the CSSI Post Action stops at the step which shows failed as status. The reboot is than not done automatically, please perform the following workaround in this case. Run the CSSI Post Actions a second time. ■ Close the CSSI Post Actions by pressing the X. ■ Start a cmd tool window (Start ⇒ All Programs⇒ Command Prompt). ■ In the cmd tool window change the directory to C:\MedCom\bin ■ Execute again the CSSI Post actions by typing in the command: <code>DiskImageRestorePostActions.exe ↵</code> Now the CSSI Post actions starts again. This takes about 20 minutes. After this time a reboot is done automatically. If the CSSI Post actions tool fails again, please contact the HSC MR in Erlangen.	

3.1.4.2 Manual Restore of System Images from USB-Disk

This procedure describes how to restore a complete software image backup or MR AWP (MRC) / MR WP (MRSC) (all hard disks and partitions) from a USB-disk.



Prerequisites



- USB-disk with actual COBR image backup
- Autologin is disabled
- Password for syadmin is set and known



Procedure if login as sysadmin is possible




- COBR image restore setps




ACTION	MRAWP/MRWP
■ Switch the MR measurement system to standby. - the syngo main UI, select System ⇒ Control. - In the System Manager pop-up, select the tab Meas & Recon. - Press the Standby button under Meas & Recon Control. - Confirm the displayed pop-up with "Yes".	☐
■ Wait until the MR measurement system is switched off.	☐
■ Remove the power line connector from the MRIR in the Electronic Power Cabinet.	☐

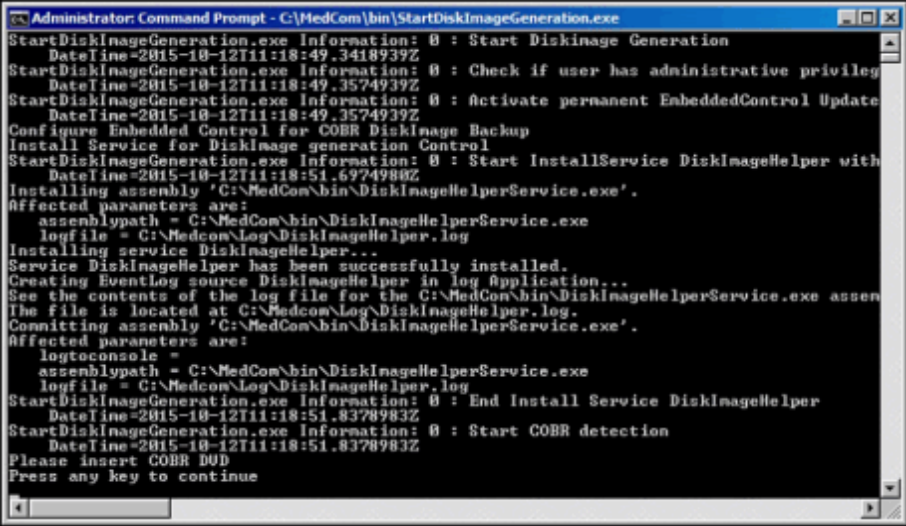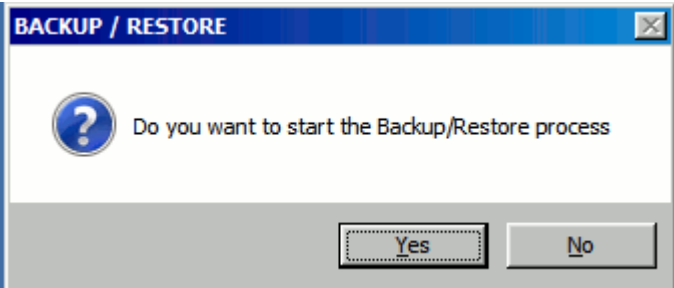
ACTION	MRAWP/MRWP
■ Reboot and log in as syadmin	☐
■ Invoke a command shell (Start ⇒ All Programs ⇒ Accessories ⇒ Command Prompt)	☐
■ In the command shell, type in: C:\MedCom\bin\StartDiskImageGeneration.exe ↵	☐
Fig. 74: Start COBR Image Creation  The command above will switch EC to permanent update mode and prepare the system for COBR image creation.	
■ The COBR tool DVD is requested. Insert the first DVD from the last created backup image into the CD/DVD-ROM drive.	☐
■ Press any key e.g. <SPACE> to start the COBR image	☐
■ Confirm the displayed popup with Yes.	☐
Fig. 75: Backup/Restore popup	
	
■ The running script performs a reboot at the end. This may take some time. Please wait!.	☐
■ A reboot is performed and the COBR tool boots.	☐

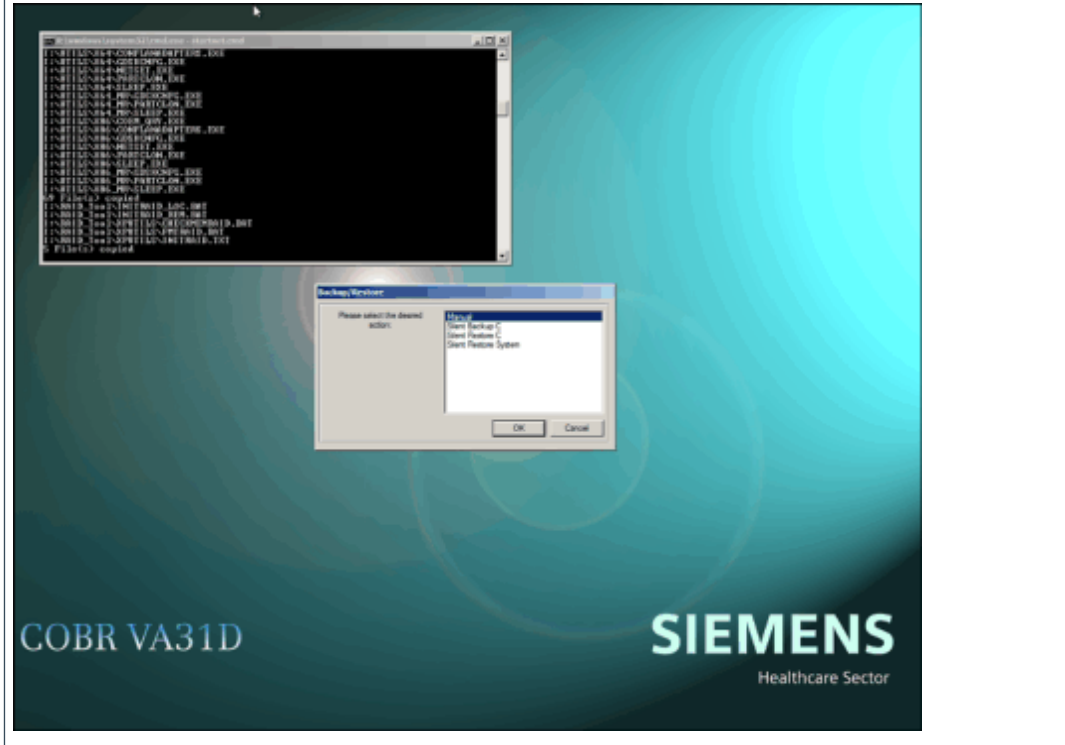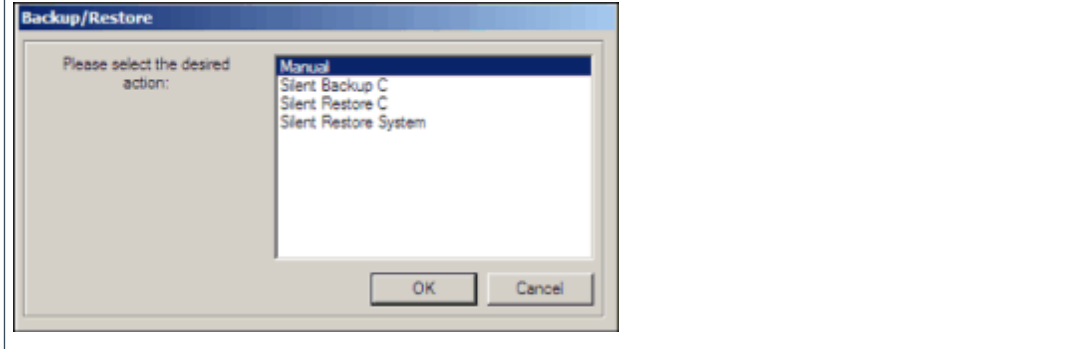ACTION	MRAWP/MRWP
Wait until the COBR tool window is displayed. Fig. 76: COBR main Backup/Restore popup 	☐
Connect the USB medium which contains the COBR image.	☐
In the Backup/Restore pop-up select <code>Manual</code> and press <code>OK</code> Fig. 77: Manual Backup/Restore popup 	☐

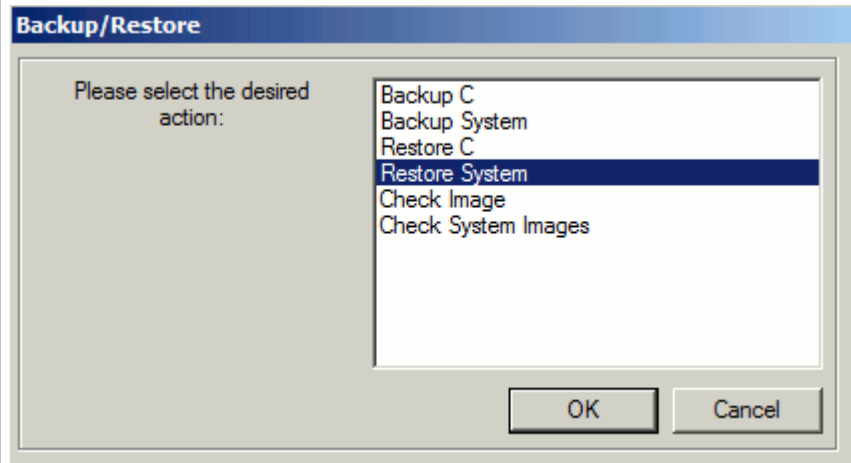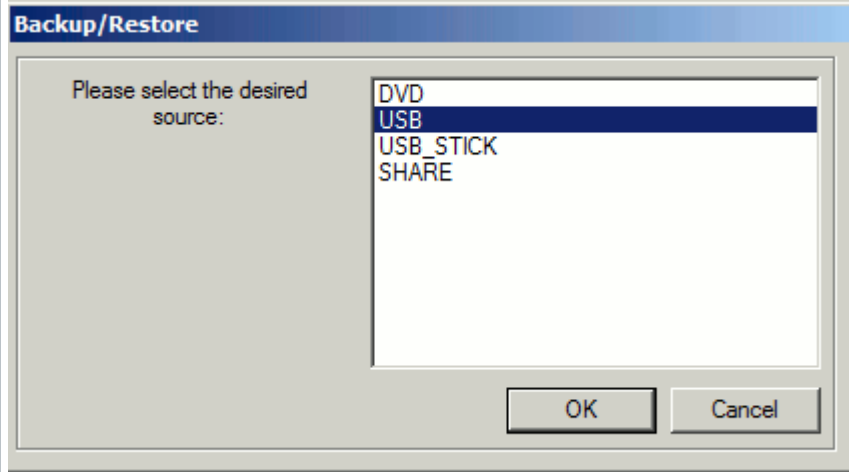ACTION	MRAWP/MRWP
■ In the next Backup/Restore pop-up select <code>Restore System</code> and press OK. <i>Fig. 78: Backup/Restore popup - Restore System</i> 	☐
■ In the next Backup/Restore pop-up select <code>USB</code> and press OK. <i>Fig. 79: Backup/Restore popup - select USB as source</i> 	☐

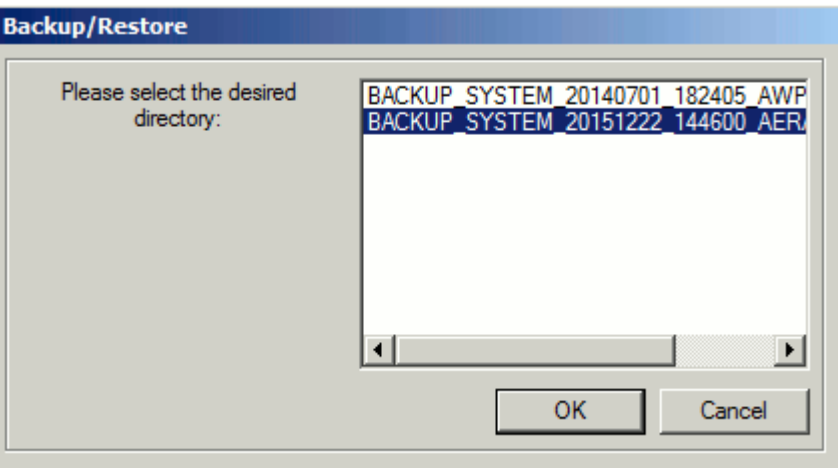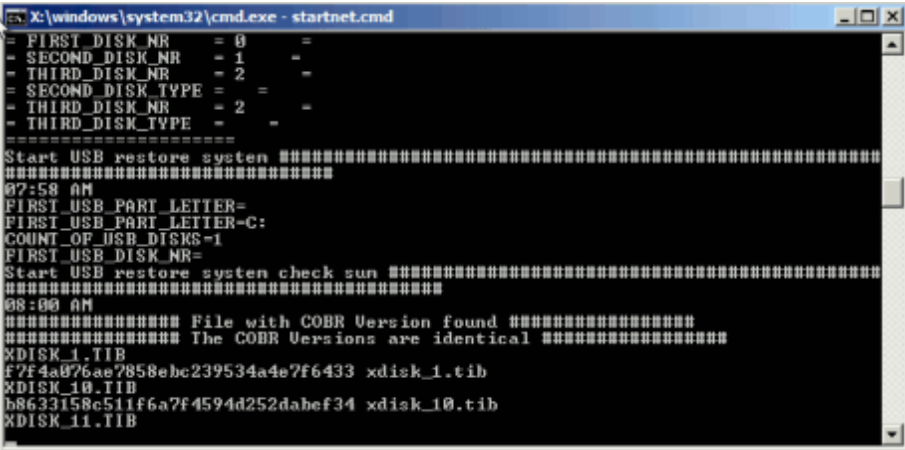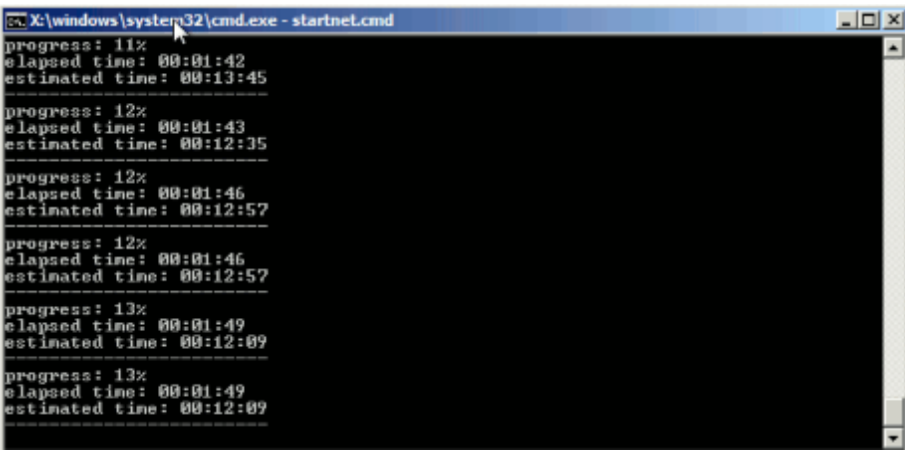ACTION	MRAWP/MRWP
In the next Backup/Restore popup select the COBR directory which you want to restore (this directory contains the desired images) and press OK. Fig. 80: Example for the image restore selection 	☐
Now the system restore will run. First the md5 verification is performedThe following is an example of the screen display. Fig. 81: COBR md5 check example output  Fig. 82: COBR image restore example output 	☐

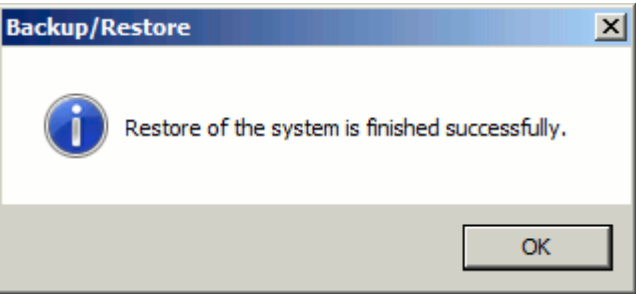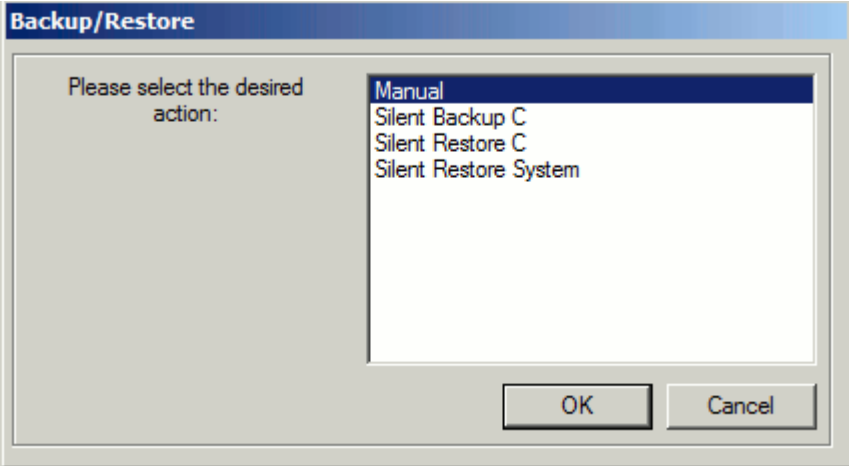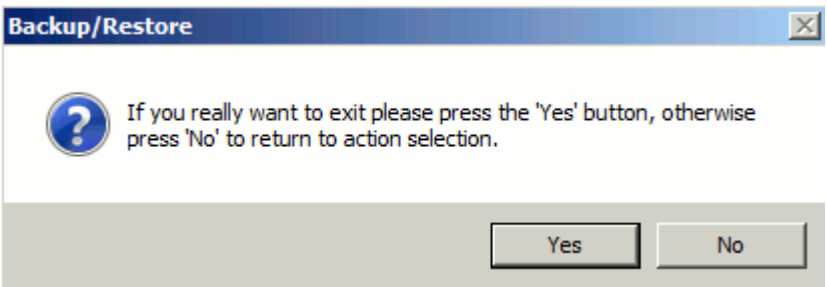ACTION	MRAWP/MRWP
Confirm the success popup by pressing OK. <i>Fig. 83: Restore finished successfully</i> 	☐
The main Backup/Restore popup will appear. <i>Fig. 84: COBR main Backup/Restore popup</i> 	☐
Press Cancel to exit the Backup/Restore procedure.	☐
Confirm the following popup with Yes. <i>Fig. 85: Confirmation popup</i> 	☐

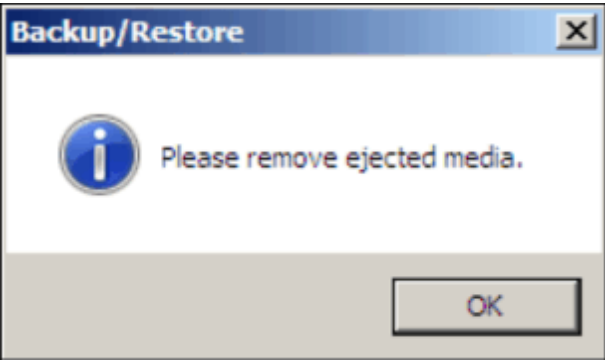ACTION	MRAWP/MRWP
The CD/DVD-ROM drive opens and the following popup is displayed. <i>Fig. 86: Backup/Restore remove ejected medium popup</i> 	☐
Remove the ejected disk from the drive.	☐
Close the drive and click OK to confirm.	☐
The system is rebooted.	☐
Disconnect the USB medium.	☐

- COBR post restore steps on MRAWP (MRC)


ACTION	MRAWP
Note: The following behaviour is new in the VB19B COBR image restore.	
After the reboot, the autologin fails with an error message. Please wait until a further reboot is done (taks about 30 seconds). The next autologin is successful done. Note: Please don't klick or hit any key ! Just wait.	☐
After the autologin, additional configuration steps are performed. Wait until they are finished, and the system is rebooted. This will take about 10 minutes. Note: Don't accept any displayed popup or press a button on any displayed popup.	☐

ACTION

MRAWP


During these steps, do not close any displayed command windows. Also do not confirm any displayed popup windows. The following windows will be displayed:

Fig. 87: Displayed popup

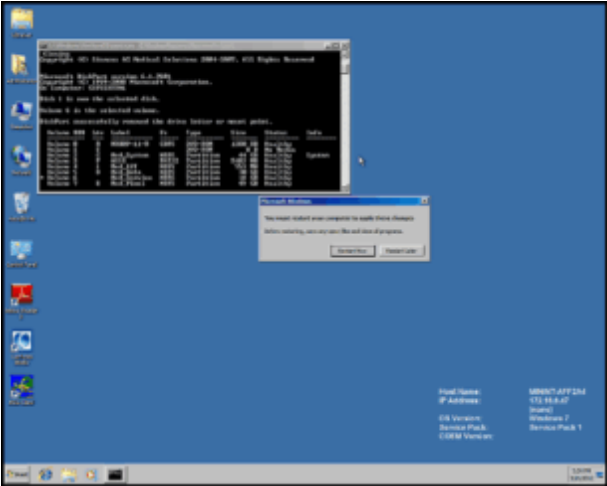


Note: Don't accept any displayed popup or press a button on any displayed popup. If the following popup is displayed, don't touch or accept it. Keep it displayed on the screen. It will be closed automatically.


Fig. 88: Windows popup

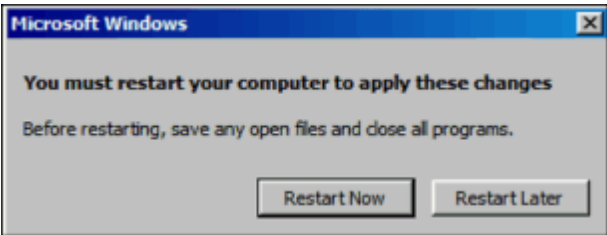

Note: The following popup can be ignored too, if it is displayed.

Fig. 89: Network location popup

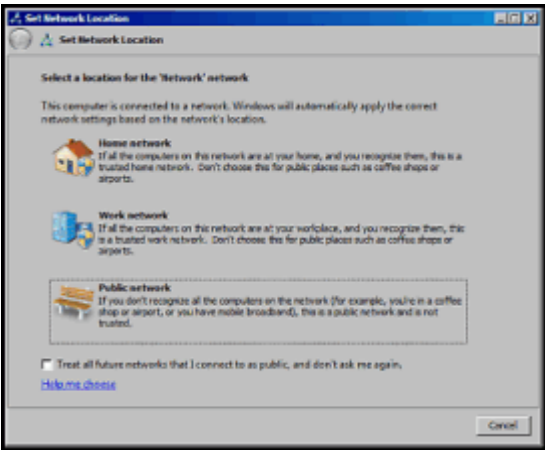

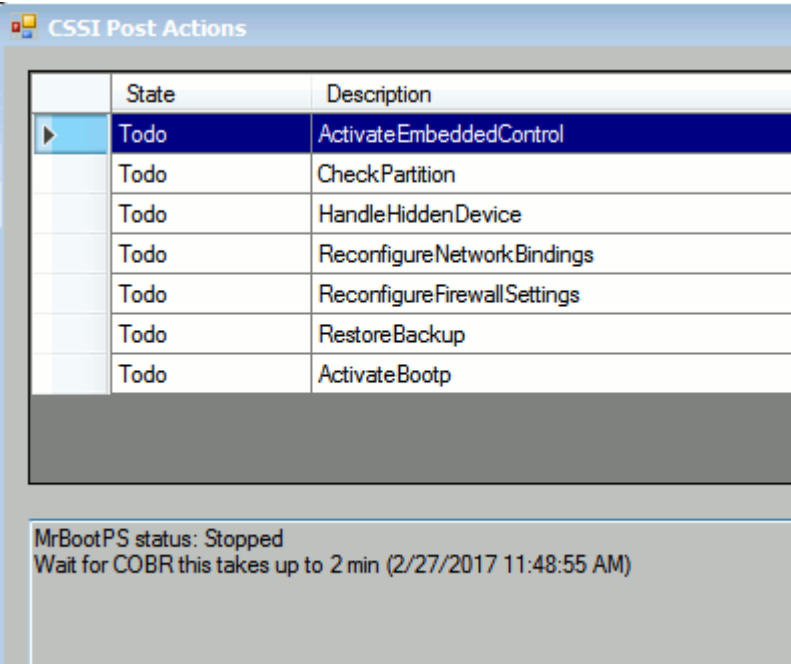ACTION	MRAWP
The CSSI Post popup is displayed. Fig. 90: CSSI Post Action popup  The following CSSI Post Action steps are performed: ■ ActivateEmbeddedControl (duration under 1 min.) ■ CheckPartition (duration under 1 min.) ■ HandleHiddenDevice (duration under 1 min.) ■ ReconfigureNetworkBindings (duration under 1 min.)	
■ If hidden network devices are found and removed, the system is rebooted automatically. After the autologin as syadmin the CSSI Post Actions popup (→ Fig. 90 Page 61) is displayed again.	☐
■ The following CSSI Post Action steps are performed: - ReconfigureFirewallSettings (duration under 1 min.) - RestoreBackup (only IP settings, duration about 1 min.) - ActivateBootp (duration under 1 min.)	☐
■ The system is rebooted automatically.	☐
■ Connect the power line connector from the MRIR in the Electronic Power Cabinet and press the SystemON button on the alarmbox.	☐
■ A autologin as meduser is done.	☐

ACTION	MRAWP
Wait until the main UI is up.	☐
Note: In case of an error the CSSI Post Action stops at the step which shows failed as status. The reboot is then not done automatically, please perform the following workaround in this case. Run the CSSI Post Actions a second time. - Close the CSSI Post Actions by pressing the X. - Start a cmd tool window (Start ⇒ All Programs ⇒ Command Prompt). - In the cmd tool window change the directory to C:\MedCom\bin - Execute again the CSSI Post actions by typing in the command: <code>DiskImageRestorePostActions.exe ↵</code> Now the CSSI Post actions starts again. This takes about 20 minutes. After this time a reboot is done automatically. If the CSSI Post actions tool fails again, please contact the HSC MR in Erlangen.	

- COBR post restore steps on MRWP (MRSC)


ACTION	MRWP
Note: The following behaviour is new in the VB19B COBR image restore.	
After the reboot, the autologin fails with an error message. Please wait until a further reboot is done (takes about 1 minute) The next autologin is successful done. Note: Please don't click or hit any key ! Just wait.	☐
After the autologin, additional configuration steps are performed. Wait until they are finished, and the system is rebooted. This will take about 5 minutes. Note: Don't close or confirm any displayed popups.	☐

ACTION

MRWP

During these steps, do not close any displayed command windows. Also do not confirm any displayed popup windows. The following windows will be displayed:


Fig. 91: Displayed popup


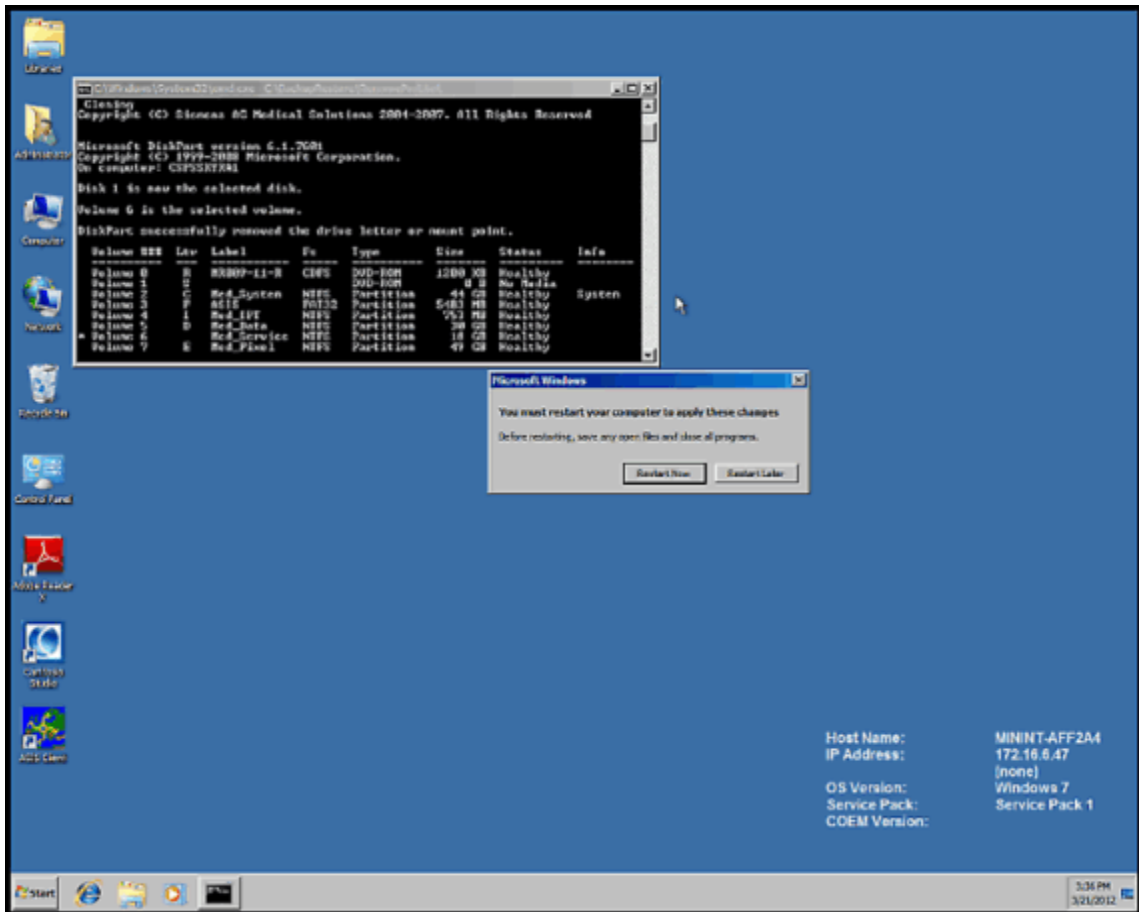


The CSSI Post popup is also displayed (see Fig. 60 “CSSI Post Actions popup”)

Note: Don’t accept any displayed popup or press a button on any displayed popup. If the following popup is displayed, don’t touch or accept it. Keep it displayed on the screen. It will be closed automatically.

Fig. 92: Windows popup

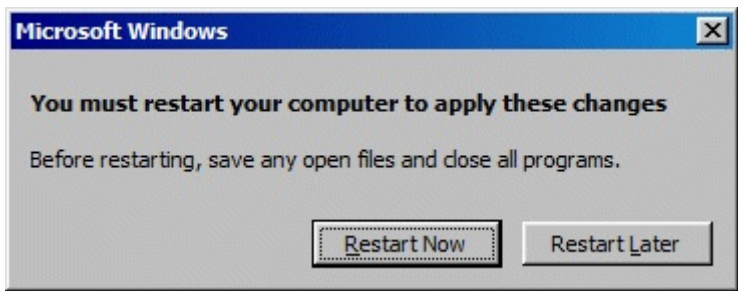

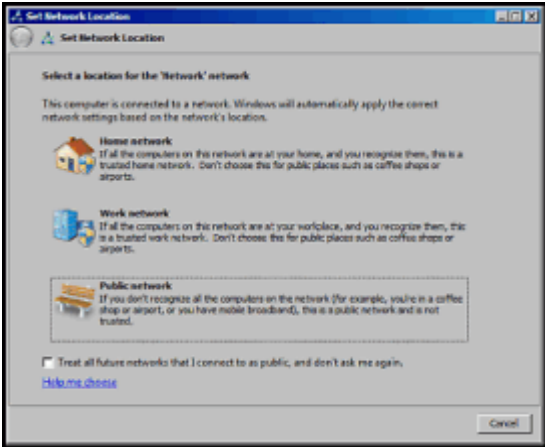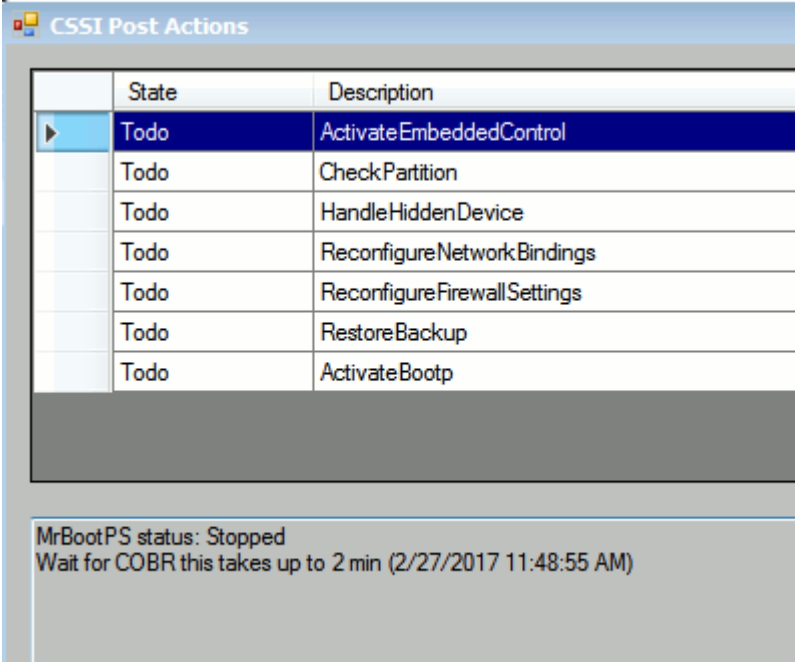ACTION		MRWP
Note: The following popup can be ignored too, if it is displayed. Fig. 93: Network location popup 		
■ The system is rebooted.		☐
■ A autologin as user "syadmin" is performed and the following CSSI post popup is displayed.		☐
■ The CSSI post popup shows the status and progress of the different post CSSI update steps. Fig. 94: CSSI Post Actions popup MRWP  Please wait until the steps are finished. During these post CSSI update steps the host will be rebooted automatically. Note: No user interaction is necessary.		☐
■ The system performs a autologin as meduser after all post CSSI update steps are done.		☐

ACTION	MRWP
Wait until the syngo main UI is up and running.	☐
Note: In case of an error the CSSI Post Action stops at the step which shows failed as status. The reboot is than not done automatically, please perform the following workaround in this case. Run the CSSI Post Actions a second time. - Close the CSSI Post Actions by pressing the X. - Start a cmd tool window (Start ⇒ All Programs⇒ Command Prompt). - In the cmd tool window change the directory to C:\MedCom\bin - Execute again the CSSI Post actions by typing in the command: <code>DiskImageRestorePostActions.exe ↵</code> Now the CSSI Post actions starts again. This takes about 20 minutes. After this time a reboot is done automatically. If the CSSI Post actions tool fails again, please contact the HSC MR in Erlangen.	



Check that the system time, time zone and date is correct.

3.1.4.3 Manual Restore of System Images from DVD

This procedure describes how to restore the complete software system for MR AWP (MRC) / MR WP (MRSC) (all hard disks and partitions) from DVD medias


Prerequisites



- Actual COBR system image on DVD is available


Procedure



- COBR image restore setps




ACTION	MRAWP/MRWP
- Switch the MR measurement system to standby. - the syngo main UI, select System ⇒ Control. - In the System Manager pop-up, select the tab Meas & Recon. - Press the Standby button under Meas & Recon Control. - Confirm the displayed pop-up with "Yes".	☐
Wait until the MR measurement system is switched off.	☐
Remove the power line connector from the MARS in the Electronic Power Cabinet. (The reason for it is that the MARS installation can run into a problem if it is not done)	☐

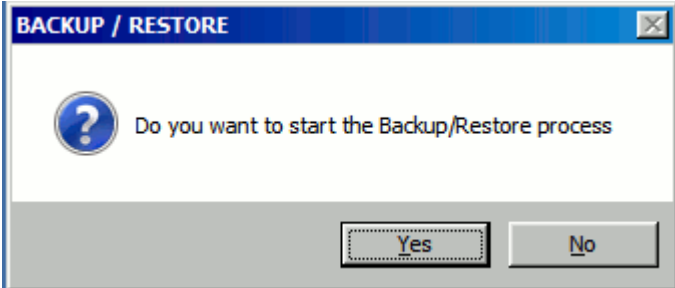
ACTION	MRAWP/MRWP
■ Reboot and log in as user administrator	☐
■ Invoke a command shell (Start ⇒ All Programs ⇒ Accessories ⇒ Command Prompt)	☐
■ In the command shell, type in: C:\MedCom\bin\StartDiskImageGeneration.exe ↵ The command above will prepare the system for COBR image restore.	☐
■ The COBR tool DVD is requested. Insert the first DVD from the last created backup image into the CD/DVD-ROM drive R:\.	☐
■ Press any key e.g.<SPACE> to start the COBR image.	☐
■ Confirm the displayed popup with Yes. <i>Fig. 95: Backup/Restore popup</i> 	☐
■ The running script performs a reboot at the end. This may take some time. Please wait!.	☐
■ A reboot is performed and the COBR tool boots.	☐

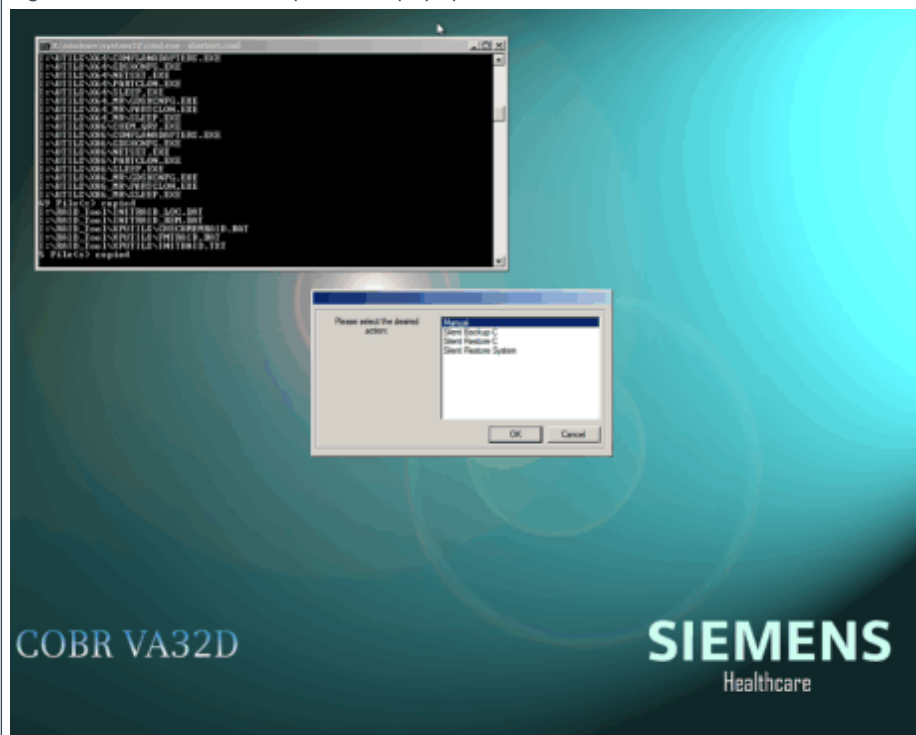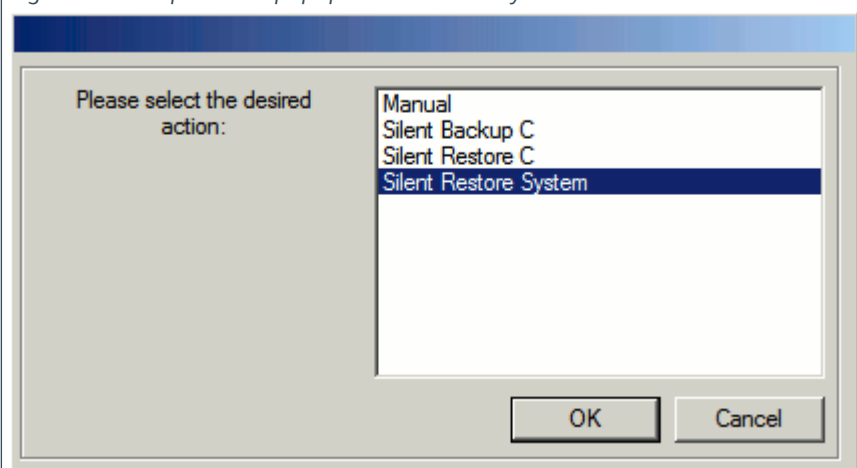ACTION	MRAWP/MRWP
■ Wait until the COBR tool window is displayed. Fig. 96: COBR main Backup/Restore popup 	☐
■ In the Backup/Restore popup, select <code>Silent Restore System</code> and press OK Fig. 97: Backup/Restore popup Silent Restore System 	☐

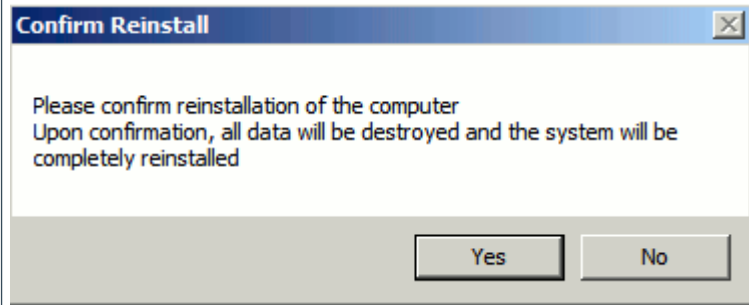ACTION	MRAWP/MRWP
Confirm the pop-up "Confirm Reinstall" with YES <i>Fig. 98: Backup/Restore popup Confirm Reinstall</i> 	☐
The system restore process starts now. The following shows an example of the screen display. <i>Fig. 99: COBR tool window</i> 	☐

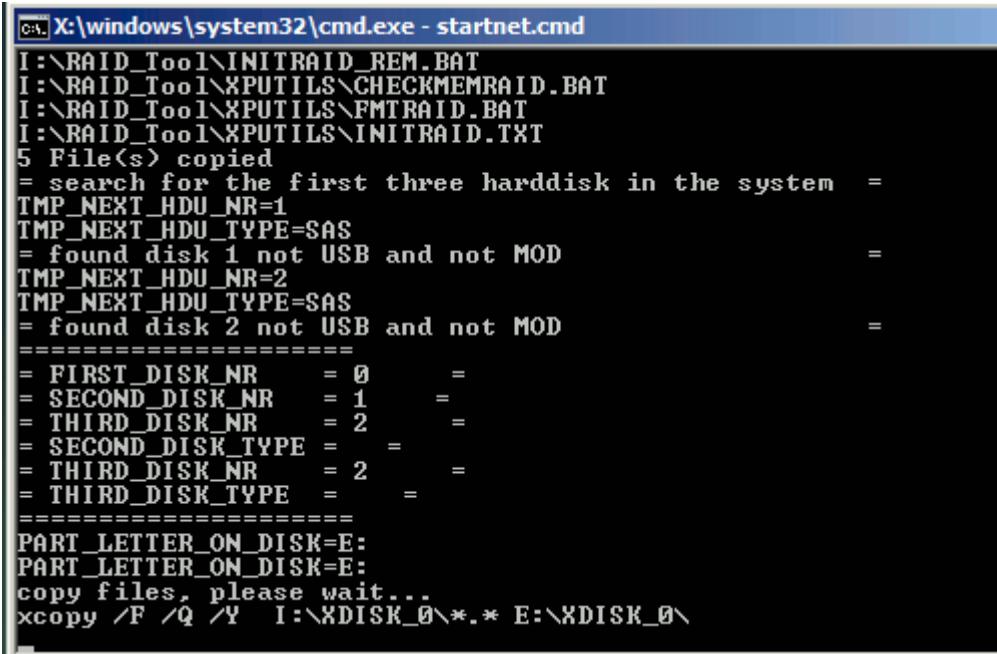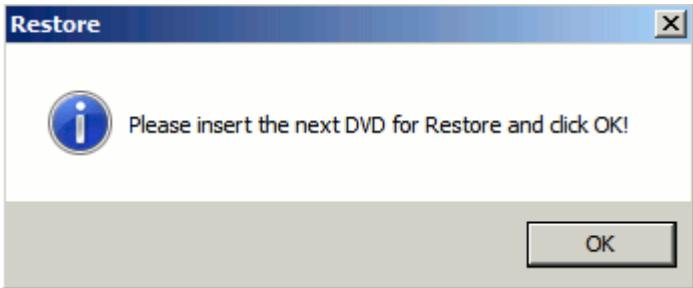ACTION	MRAWP/MRWP
First all image files from the DVD for the system disk (disk1) are copied from the DVD to the hard disk. Fig. 100: Copy image files <pre>C:\X:\windows\system32\cmd.exe - startnet.cmd I:\RAID_Tool\INITRAID_REM.BAT I:\RAID_Tool\XPUTILS\CHECKMEMRAID.BAT I:\RAID_Tool\XPUTILS\FMTRAID.BAT I:\RAID_Tool\XPUTILS\INITRAID.TXT 5 File(s) copied = search for the first three harddisk in the system = TMP_NEXT_HDU_NR=1 TMP_NEXT_HDU_TYPE=SAS = found disk 1 not USB and not MOD = TMP_NEXT_HDU_NR=2 TMP_NEXT_HDU_TYPE=SAS = found disk 2 not USB and not MOD = ===== = FIRST_DISK_NR = 0 = = SECOND_DISK_NR = 1 = = THIRD_DISK_NR = 2 = = SECOND_DISK_TYPE = = = THIRD_DISK_NR = 2 = = THIRD_DISK_TYPE = = ===== PART_LETTER_ON_DISK=E: PART_LETTER_ON_DISK=E: copy files, please wait... xcopy /F /Q /Y I:\XDISK_0*.* E:\XDISK_0\</pre> After the files are copied, the DVD is automatically ejected.	☐
Remove the ejected medium from the CD/DVD-ROM drive, insert the next DVD from the system image DVD set, and close the drive.	☐
Confirm the request popup with OK. Fig. 101: Next DVD request popup 	☐
The copy procedure continues. The next DVD is requested automatically until all required data are copied. The copy and the confirmation of the popup (→ Fig. 101 Page 69) have to be repeated for all available DVDs.	

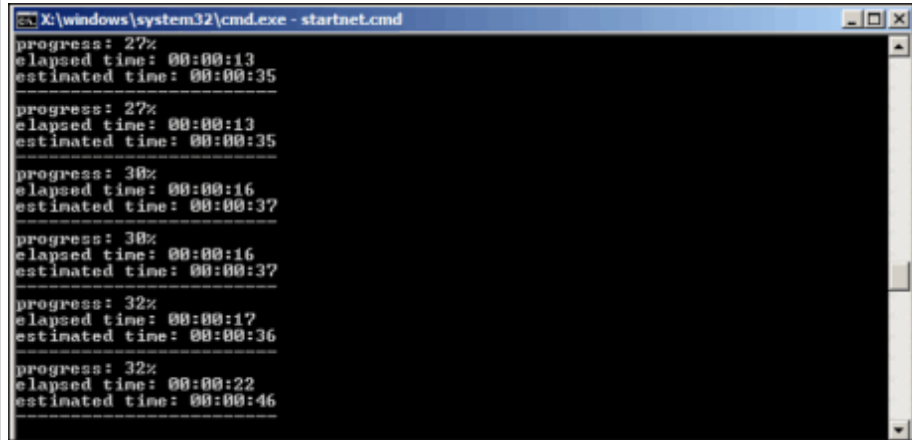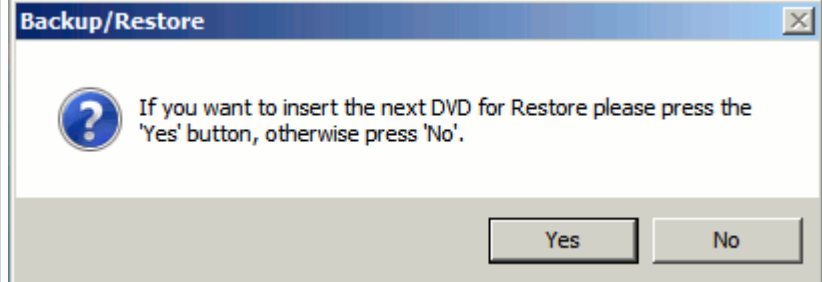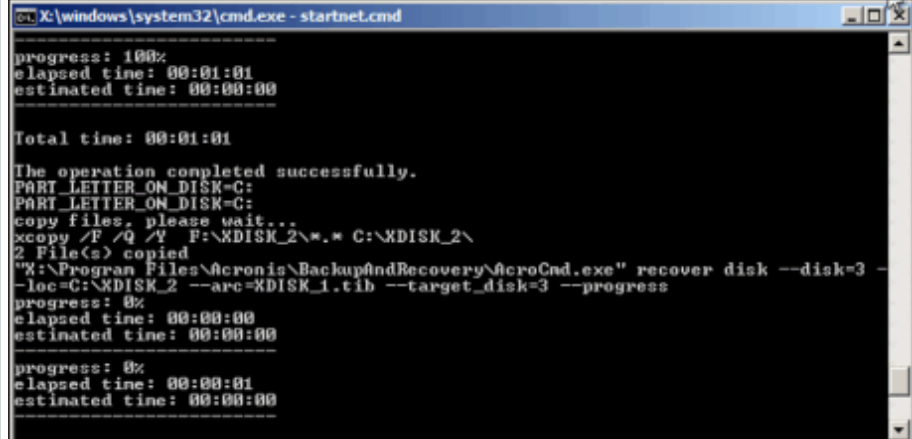ACTION	MRAWP/MRWP
After the last DVD has been copied, the Image Restore starts <i>Fig. 102: Image restore progress output</i> 	☐
Wait until the following popup appears. <i>Fig. 103: Next DVD request popup</i> 	☐
Confirm the above shown "Next DVD request popup" with "No".	☐
Now the last disk is restored from the image. <i>Fig. 104: Last disk restore</i> 	☐

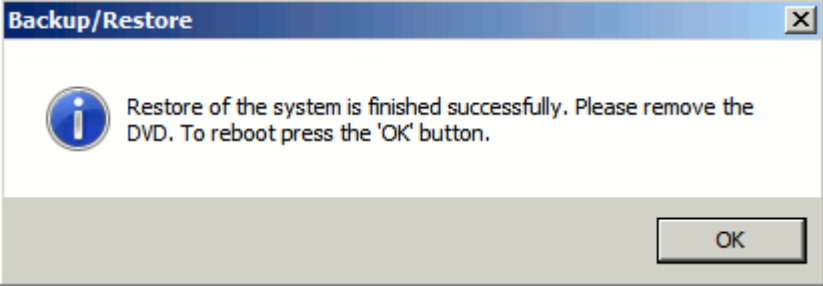ACTION	MRAWP/MRWP
When all data have been successfully restored, the following popup is displayed. <i>Fig. 105: Backup/Restore successfully finished popup</i> 	☐
Remove the ejected medium, close the drive, and confirm this popup with OK.	☐
The system is rebooted.	☐

- COBR post restore steps on MRAWP (MRC)


ACTION	MRAWP
After the autologin, additional configuration steps are performed. Wait until they are finished, and the system is rebooted. This will take about 10 minutes. Note: Don't accept any displayed popup or press a button on any displayed popup.	☐

ACTION

MRAWP


During these steps, do not close any displayed command windows. Also do not confirm any displayed popup windows. The following windows will be displayed:

Fig. 106: Displayed popup

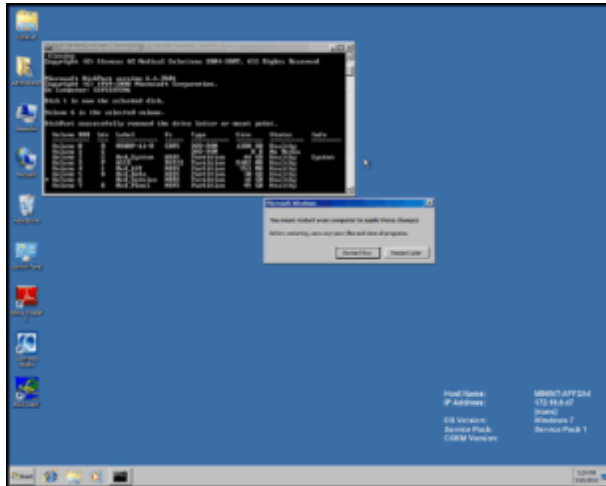


Note: Don't accept any displayed popup or press a button on any displayed popup. If the following popup is displayed, don't touch or accept it. Keep it displayed on the screen. It will be closed automatically.

Fig. 107: Windows popup

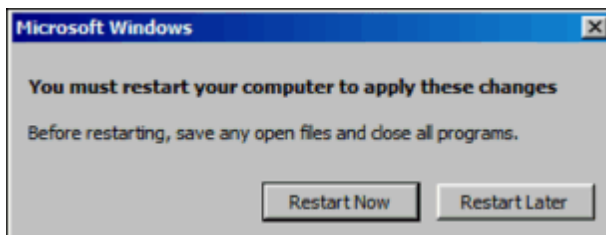

Note: The following popup can be ignored too, if it is displayed.

Fig. 108: Network location popup

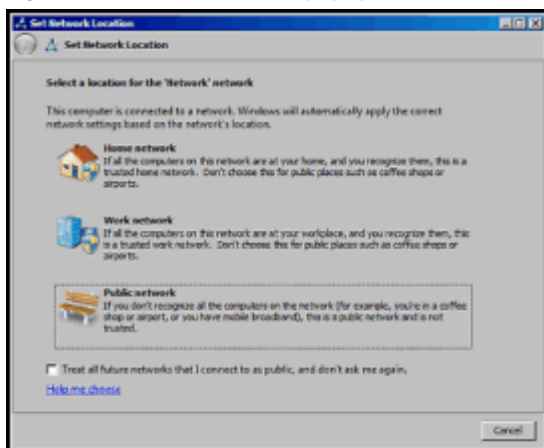


ACTION

MRAWP


The CSSI Post popup is displayed.

Fig. 109: CSSI Post Action popup

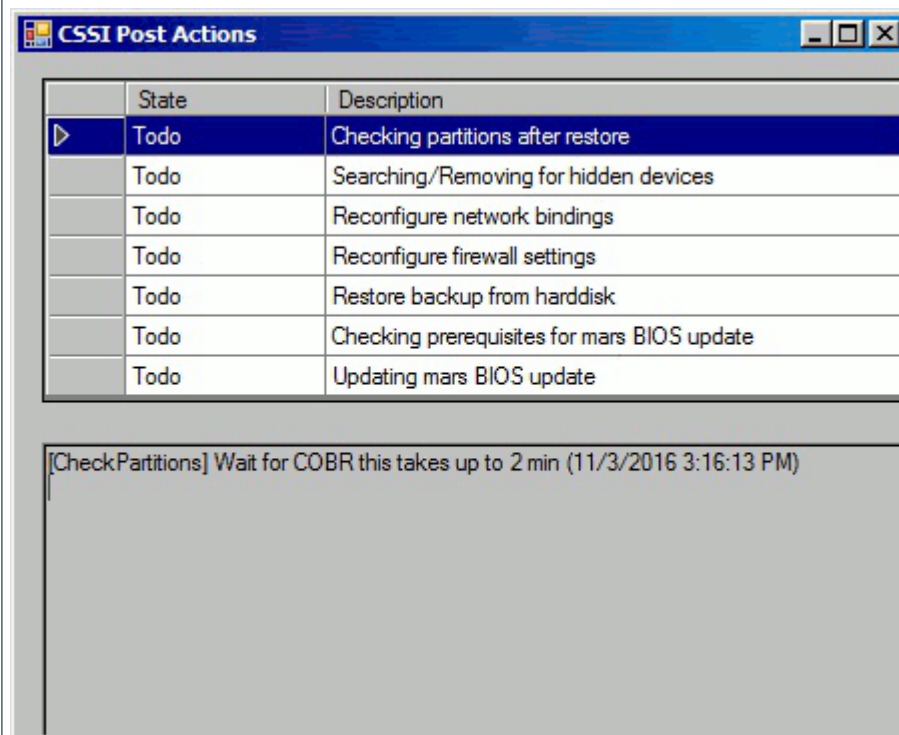

The following CSSI Post Action steps are performed:

- Checking partitiony after restore (duration under 1 min.)
- Searching/removeing hidden device (duration under 1 min.)
- Reconfigure network bindings (duration under 1 min.)
- If hidden network devices are found and removed, the system is rebooted automatically.
After the autologin as **administrator** the CSSI Post Actions popup (→ Fig. 109 Page 73) is displayed again. ☐
- The following CSSI Post Action steps are performed:
 - Reconfigure firewall settings (duration under 1 min.)
 - Restore backup (only IP settings, duration about 1 min.)☐
- The system is rebooted automatically. ☐

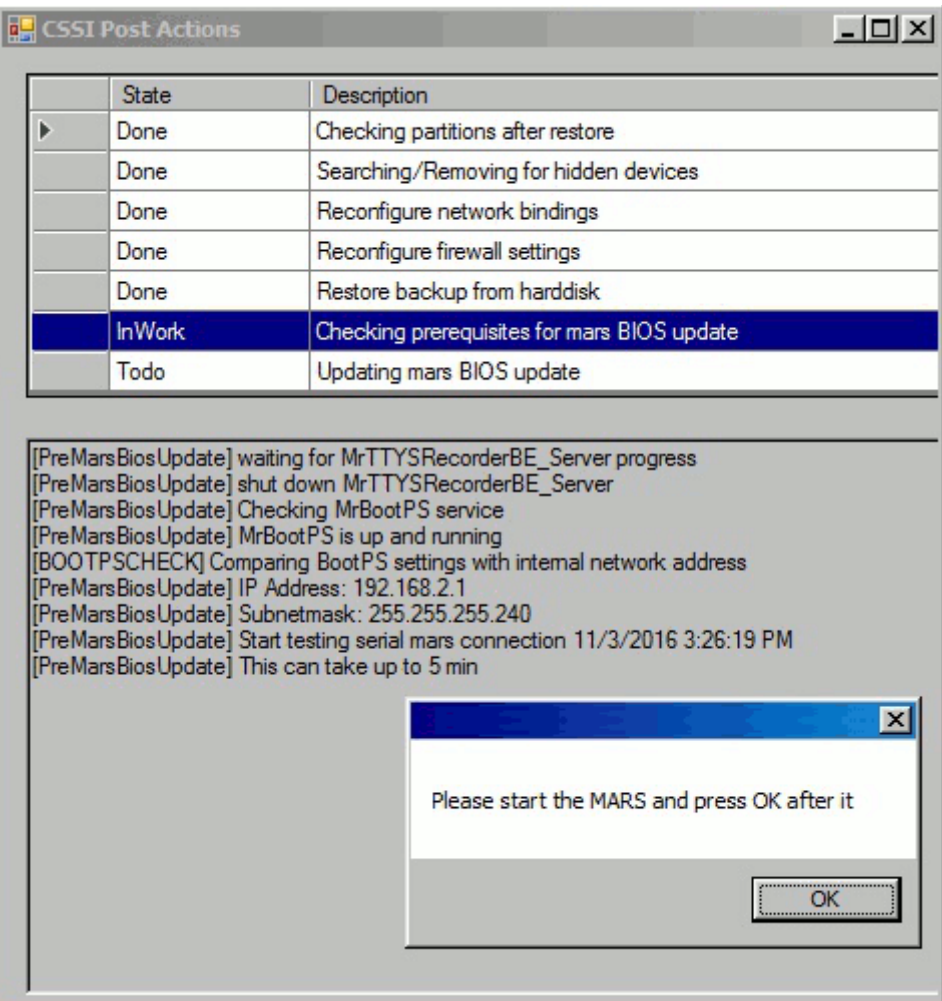ACTION	MRAWP
■ A autologin as user "administrator" is performed and the CSSI Post Actions popup (→ Fig. 109 Page 73) is displayed. The following CSSI Post Action steps are performed: - Install Mars Wait until you get the Please start the MARS... popup (duration about 5 min.). <i>Fig. 110: MARS Start popup</i> 	☐
■ Connect the power line connector from the MARS in the Electronic Power Cabinet and press the SystemON button on the alarmbox.	☐
■ Confirm the Please start the MARS... popup with OK.	☐
■ The MARS is booted up and the following CSSI post action step continues: - Install Mars (duration about 15 min.) Automatically a reboot is performed after this step is done.	☐


ACTION	MRWP
Note: In case of an error the CSSI Post Action stops at the step which shows failed as status. The reboot is than not done automatically, please perform the following workaround in this case. Run the CSSI Post Actions a second time. ■ Close the CSSI Post Actions by pressing the X. ■ Start a cmd tool window (Start ⇒ All Programs⇒ Command Prompt). ■ In the cmd tool window change the directory to C : \MedCom\bin ■ Execute again the CSSI Post actions by typing in the command: <code>CSSIPostRestoreProcess.exe ↵</code> Now the CSSI Post actions starts again. This takes about 20 minutes. After this time a reboot is done automatically. If the CSSI Post actions tool fails again, please contact the HSC MR in Erlangen.	
■ Wait until the main UI is up.	☐

- COBR post restore steps on MRWP (MRSC)

ACTION	MRWP
■ After the autologin, additional configuration steps are performed. Wait until they are finished, and the system is rebooted. This will take about 5 minutes. Note: Don't close or confirm any displayed popups.	☐

During these steps, do not close any displayed command windows. Also do not confirm any displayed popup windows. The following windows will be displayed:


Fig. 111: Displayed popup


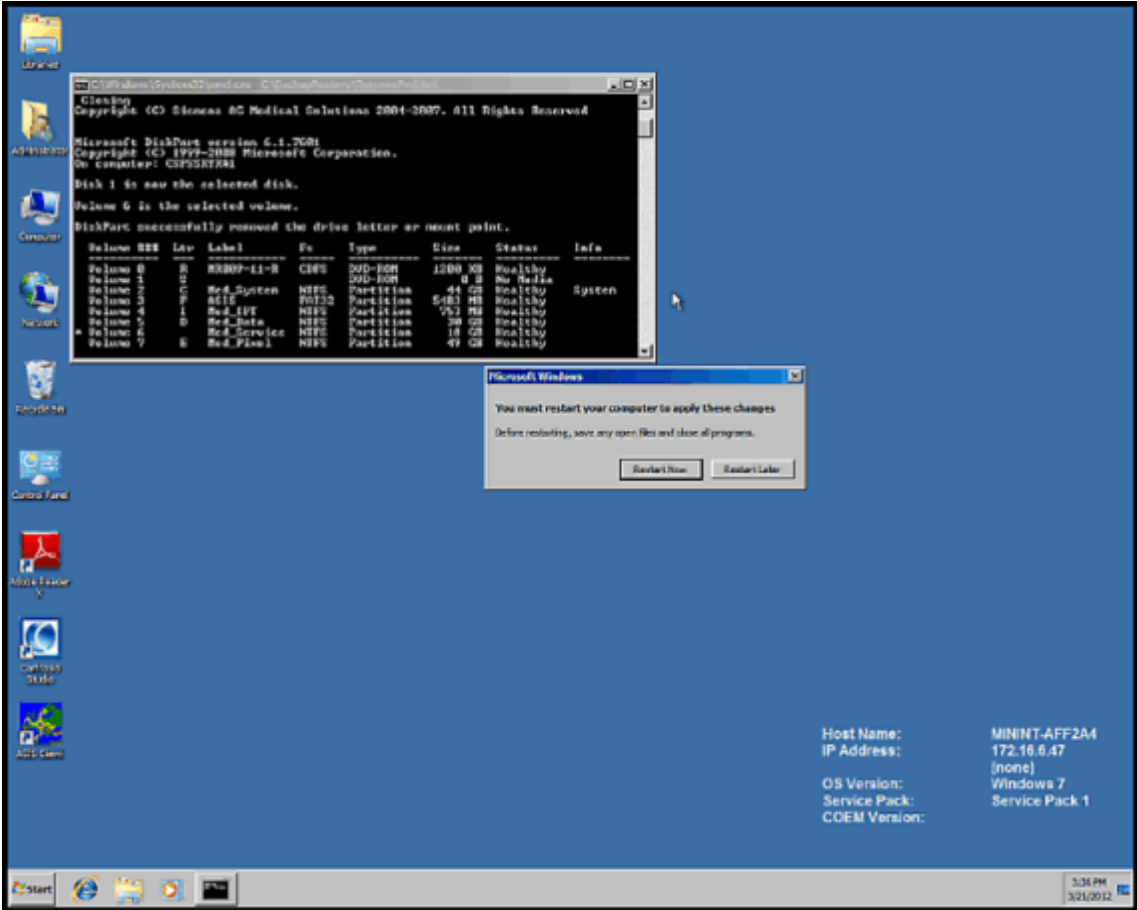


The CSSI Post popup is also displayed (see Fig. 60 “CSSI Post Actions popup”)

Note: Don’t accept any displayed popup or press a button on any displayed popup. If the following popup is displayed, don’t touch or accept it. Keep it displayed on the screen. It will be closed automatically.

Fig. 112: Windows popup

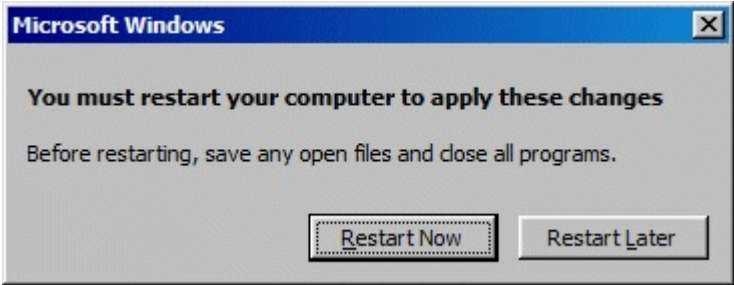

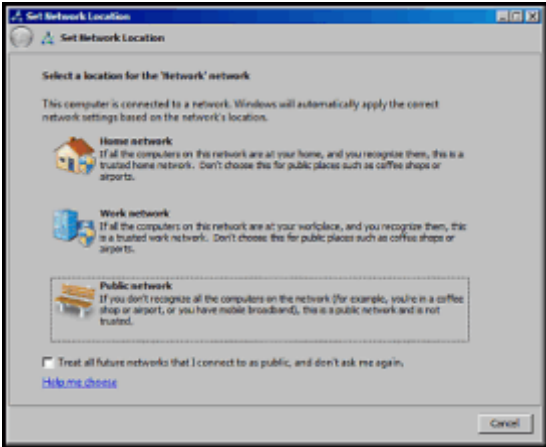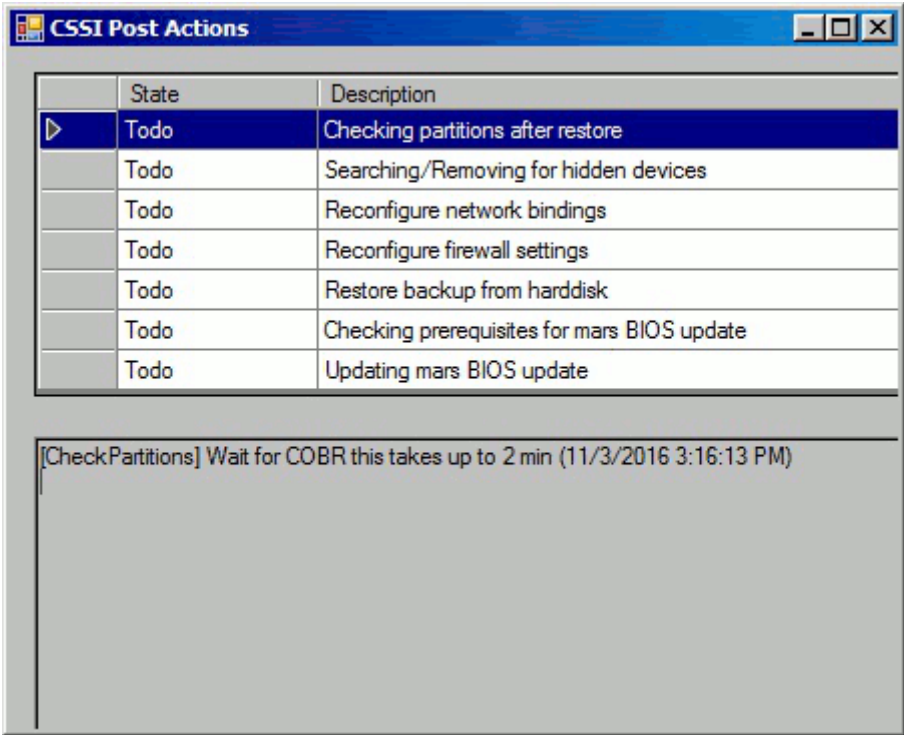ACTION		MRWP
Note: The following popup can be ignored too, if it is displayed. Fig. 113: Network location popup 		
■ The system is rebooted.		☐
■ A autologin as user "administrator" is performed and the following CSSI post popup is displayed.		☐
■ The CSSI post popup shows the status and progress of the different post CSSI up-date steps. Fig. 114: CSSI Post Actions popup 		☐
Please wait until the steps are finished. During these post CSSI update steps the host will be rebooted automatically. Note: No user interaction is necessary.		

ACTION	MRWP
■ The system performs a autologin as meduser after all post CSSI update steps are done.	☐
■ Wait until the syngo main UI is up and running.	☐



Check that the system time, time zone and date is correct.

3.1.5 Restore site-specific data from CD

Restoring the site-specific data may be helpful in many situations. Please read the following information before beginning with the work.

In case of corrupted site-specific data, we recommend saving the actual status of the site-specific data to another medium before restoring an old save set. This will save time if the problem is resolved by reinstalling the software. The actual tune-up data may save a great deal of work, especially if the customer optimized his protocols recently.

- Login as **administrator** ↵.
- Start the "Internet Explorer"
- Login to the service platform, select **Backup & Restore**.
- Insert the backup disk in the CD/DVD-ROM.
- Select **Restore** in the "Command" pull-down menu.
- Select the source drive CD/DVD-ROM (drive letter **R:**) in the "Drives" menu.
- Select your xxx_**Master-Backup**_xxxx_xxxx_ok.ar to be restored in the "Archive" menu.
- In the menu "Groups", select the necessary packages to be restored.
- Click **Go** to start the "Restore", check for errors in the log window after restore has been completed. (Estimated time ~10-15 min.)



To ensure the Administrator account password is correct following restoration of the dynamic data, it is recommended to set the default Administrator password in the Service GUI > **Configuration** > **Local Host** > **Users**



To ensure the customer username e.g. **meduser** is correct after restoration of the dynamic data, it is recommended check with the customer in the Service GUI > **Configuration** > **Local Host** > **Users**

3.1.5.1 Final checks



Check what is most recent the system "Image" or your set of "site specific" data

- Log in as **meduser** ↵.
- Open the service platform; select **Configuration > Meas. Settings** > click **Save**.
- Perform the "Quality Assurance" procedure: "Image Orientation" in the service platform..
- Make sure that the system is generating images.

Setting that are not part of the site-specific data:

Remote access mode (No Access, Limited Access, etc. in the customer UI).

Auto-transfer flag ("Configuration", "Service", "Auto Transfers")

Time settings ("Configuration", "Service", "Site Info, setp 2")

The mark on the checkbox DaylightTime automatic must be removed.

Postscript Printer Configuration (Customer paper printer)

ImageText Settings (ImageText Settings from ImageText Dialog)

Time zone

Transfer Rules

Autologin for meduser (frontend username)

Frontend username (e.g. **meduser**)

Secure DICOM printer setting (must be manually configured)

3.1.6 Restore site-specific from external USB device



Please check the external media (e.g. USB-disk) first with the virus scanner of the service PC. Be sure that the external media used is free of viruses.

Restoring the site-specific data may be helpful in many situations. Please read the following information before beginning with the work.

In case of corrupted site-specific data, we recommend saving the actual status of the site-specific data to another medium before restoring an old save set. This will save time if the problem is resolved by reinstalling the software. The actual tune-up data may save a great deal of work, especially if the customer optimized his protocols recently.

- Login as **administrator** ↵.
- Start the "Internet Explorer"
- Login to the service platform, select **Backup & Restore**.

- Connect the USB device.
- Select **Restore** in the "Command" drop-down menu.
- Select the source USB-drive (drive letter e.g. **F:**) in the "Drives" menu.
- Select your xxx_**Master-Backup**_xxxx_xxxx_ok.ar to be restored in the "Archive" menu.
- In the menu "Groups", select the necessary packages to be restored.
- Click **Go** to start the "Restore", check for errors in the log window after restore has been completed. (Estimated time ~20-60 min.)



To ensure the Administrator account password is correct following restoration of the dynamic data, it is recommended to set the default Administrator password in the Service GUI > **Configuration** > **Local Host** > **Users**



To ensure the customer username e.g. **meduser** is correct after restoration of the dynamic data, it is recommended check with the customer in the Service GUI > **Configuration** > **Local Host** > **Users**

3.1.6.1 Final checks



Check what is most recent the system "Image" or your set of "site specific" data

- Log in as **meduser** ↵.
- Open the service platform; select **Configuration** > **Meas. Settings** > click **Save**.
- Perform the "Quality Assurance" procedure: "Image Orientation" in the service platform..
- Make sure that the system is generating images.

Setting that are not part of the site-specific data:

Remote access mode (No Access, Limited Access, etc. in the customer UI).

Auto-transfer flag ("Configuration", "Service", "Auto Transfers")

Time settings ("Configuration", "Service", "Site Info, setp 2")

The mark on the checkbox DaylightTime automatic must be removed.

Postscript Printer Configuration (Customer paper printer)

ImageText Settings (ImageText Settings from ImageText Dialog)

Time zone

Transfer Rules

Autologin for meduser (frontend username)

Setting that are not part of the site-specific data:Frontend username (e.g. **meduser**)

Secure DICOM printer setting (must be manually configured)

4.1 Configuration

4.1.1 General

For detailed information regarding the configuration of syngo MR systems, refer to **Help** in the Service Software Platform

If you have changed the configuration within the Service Software Platform and the message “reboot system pending” is displayed in the “status bar”, reboot the system as follows: Click **Home**, confirm the reboot message window with **OK**.



Language settings and the configuration UI for option Chorus are described in chapter Appendix.



How to configure Siemens Remote Service is described separately in the document “Installation and Startup for SRS”.



The use of non-released combinations of hardcopy cameras and BL products is prohibited and does not meet the requirements defined. For detailed information regarding the CAMERA configuration, refer to technical documentation “Camera Info ...” valid for your NUMARIS/4 software version.



The Interoperability Competence Center provides a powerful configuration management tool, refer to:

https://intranet.healthcare.siemens.com/en/home/Pages/Healthcare_Global_Intranet.aspx

Then navigate to **DI - Diagnostic Imaging > SY > Healthcare Topics > Interoperability Competence Center**

https://intranet.healthcare.siemens.com/cms/SY/en/Med-wideTopics/Connectivity_Competence_Center/Pages/Interoperability_Competence_Center.aspx

Afterwards select **Healthcare Interoperability Database**.



Remember to create a backup if the configuration has changed.



The Normalization flag is reset after configuration changes. Please perform the Normalization after changing the configuration or restoring dynamic data.

The Normalization flag is also reset after an update of the license. Please perform the Normalization after the update.

4.1.2 Network Configuration Restrictions



In case the measurement network conflicts with the customer network, see chapter Appendix, Troubleshooting, Customer Network IP Conflicts with Measurement Network .

4.1.3 Activation of Security Settings



HIPAA configuration is only to be performed by personal trained and certified by the "Siemens Medical Academy".

During the first configuration of the security settings following error message is displayed: "retrieve of security data failed". After confirming the pop-up with OK the checkboxes in the security properties will be dimmed.

Workaround:

- Start the local service UI at the MRC
- Reboot the system
- Configure the **Security Settings**
- Re-boot the system once more for activation of the security settings

5.1 Appendix

5.1.1 Monitor Settings

Check the monitor settings in the local service menu.


- Select: Options -> Local Service

or as user administrator (Start the Internet Explorer and type in the following URL to start the local service UI: <http://localhost/service>)

login (service key is needed) select "Configuration" check the point "MonitorType"; if you don't have the option "Multi Monitor" (this is not the option syngo.via) the checkbox "Syngo single Monitor" has to be selected!


If you have the option syngo.via the the checkbox "Syngo single Monitor" has to be selected too!.

Fig. 115: Monitor settings

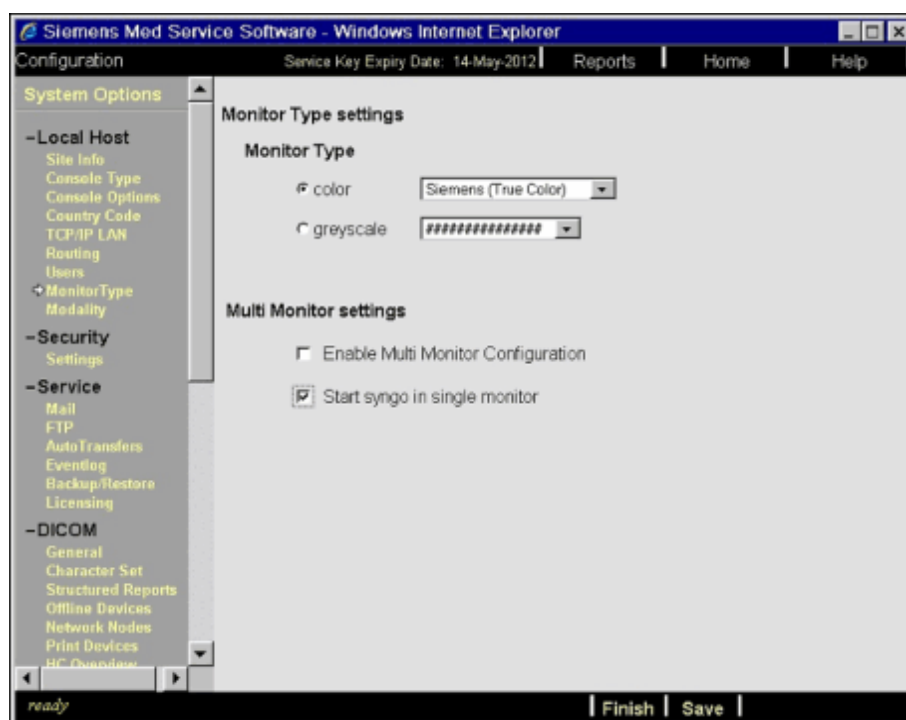

- If the display order or ScreenResolution is incorrect, select in the Control Panel 'Display' 'ScreenResolution' and set the correct order/ScreenResolution.



Fig. 116: Display resolution

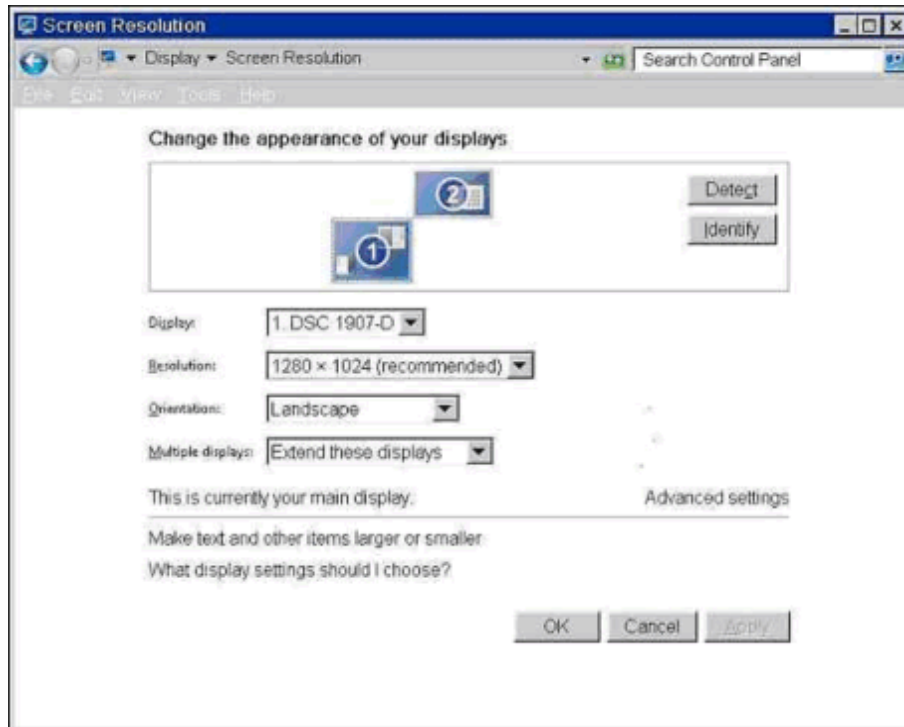

- Main display resolution: 1280 x 1024
- Patient Data Display (PDD) resolution : 1280 x 800

5.1.2 License update

ACTION	MRC	MRSC
■ Be sure to be logged in as meduser.	☐	☐
■ Go to local service.	☐	☐
■ Select Configuration in the service platform.	☐	☐
■ Insert the media with the new license file	☐	☐
■ Select Licensing and import the new license file.	☐	☐

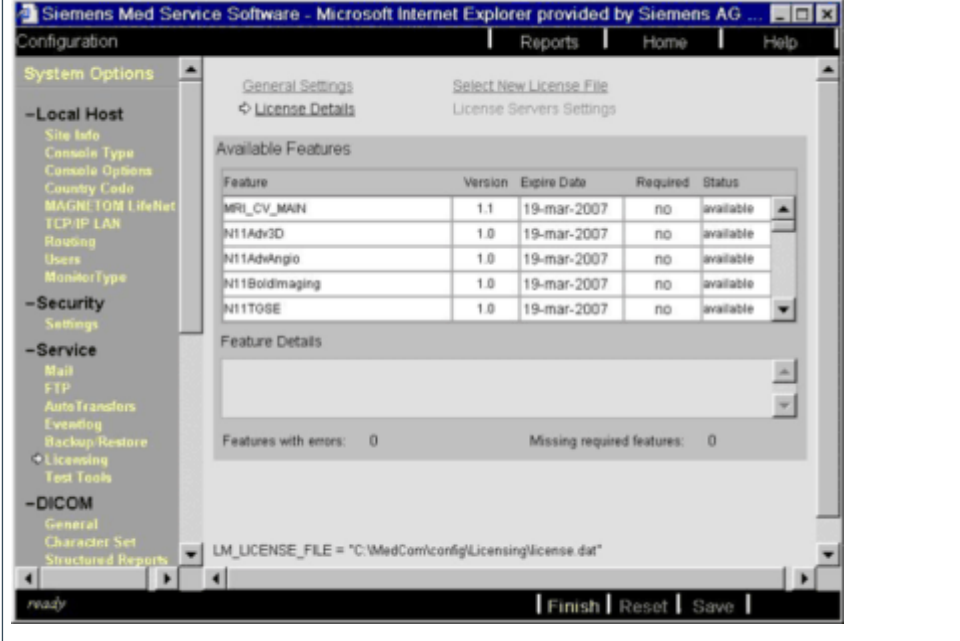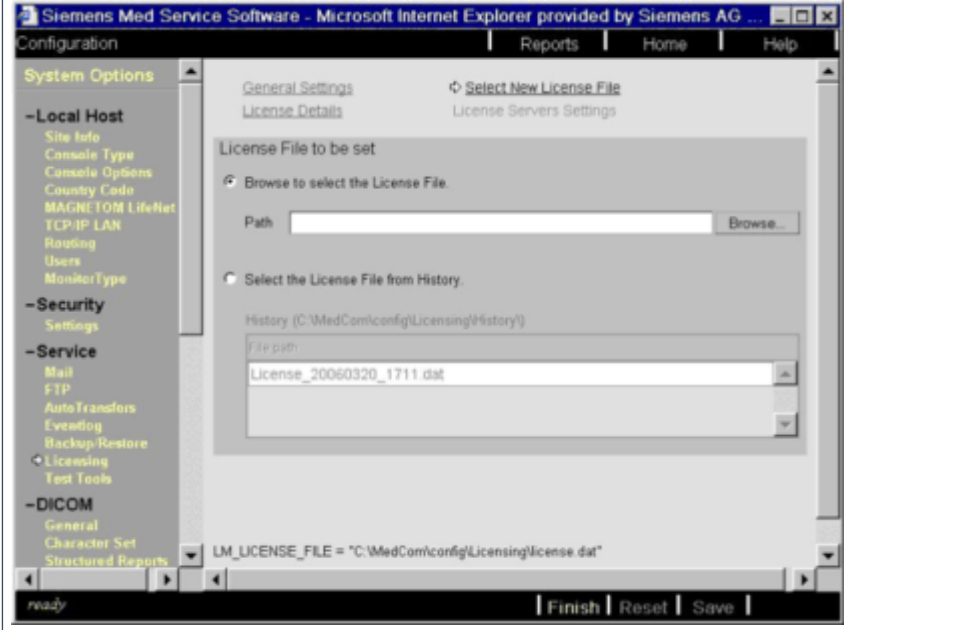ACTION	MRC	MRSC
■ In the Licensing UI, choose "Select New License File" Fig. 117: The screenshot shows the 'Select New License File' tab in the Siemens Med Service Software configuration utility. The left sidebar contains a tree view with categories like Local Host, Security, Service, and DICOM. The main area has tabs for General Settings, License Details, and Select New License File (which is active). Under 'Select New License File', there are sub-tabs for License Servers Settings and Available Features. The 'Available Features' section contains a table with columns: Feature, Version, Expire Date, Required, and Status. The table lists several features like MRI_CV_MAIN, N11Adv3D, N11AdvAngio, N11Boldmaging, and N11T0SE, all with version 1.0 and an expire date of 19-mar-2007. Below the table is a 'Feature Details' section and a status line indicating 'Features with errors: 0' and 'Missing required features: 0'. At the bottom, there's a field for 'LM_LICENSE_FILE' and buttons for 'Finish', 'Reset', and 'Save'.	☐	☐
■ Under "Select New License File": Fig. 118: The screenshot shows the 'Select New License File' sub-tab. It has two radio buttons: 'Browse to select the License File' (which is selected) and 'Select the License File from History'. Below the first radio button is a 'Path' text field and a 'Browse...' button. Below the second radio button is a 'History' section with a list of file paths, including 'License_20060320_1711.dat'. At the bottom, there's a field for 'LM_LICENSE_FILE' and buttons for 'Finish', 'Reset', and 'Save'. - Use the Browse button to navigate to the source for the new li- cense file - Import the new license file by clicking Save.	☐	☐


5.1.3 Restoring Customer Protocols from Last Backup

- Restoring the backup
 - After completing installation of the MRease Software, select **Options** from the menu, then select **Local Service**.
 - Enter the valid service key and click OK to confirm.
 - Insert the last created backup CD which includes the package "CustomerProtocols" into the CD/DVD writer.
 - Select **Backup/Restore**.
 - Use the **Restore** command to select the CD/DVD writer (S:I) and select the **CustomerProtocols** package.
 - Click **GO** to restore the programs.
 - Click **Home** and perform the requested reboot.


5.1.4 Dongle exchange and license update

ACTION	MRC	MRSC
■ Shut down the system and switch it off	☐	n.a.
■ Shut down the MRSC and switch it off	n.a.	☐
■ Replace the old dongle with the new one.	☐	☐
■ Switch on the system.	☐	☐
■ Be sure to be logged in as meduser. If so, the following message is displayed in the IE window.	☐	☐

Fig. 119:

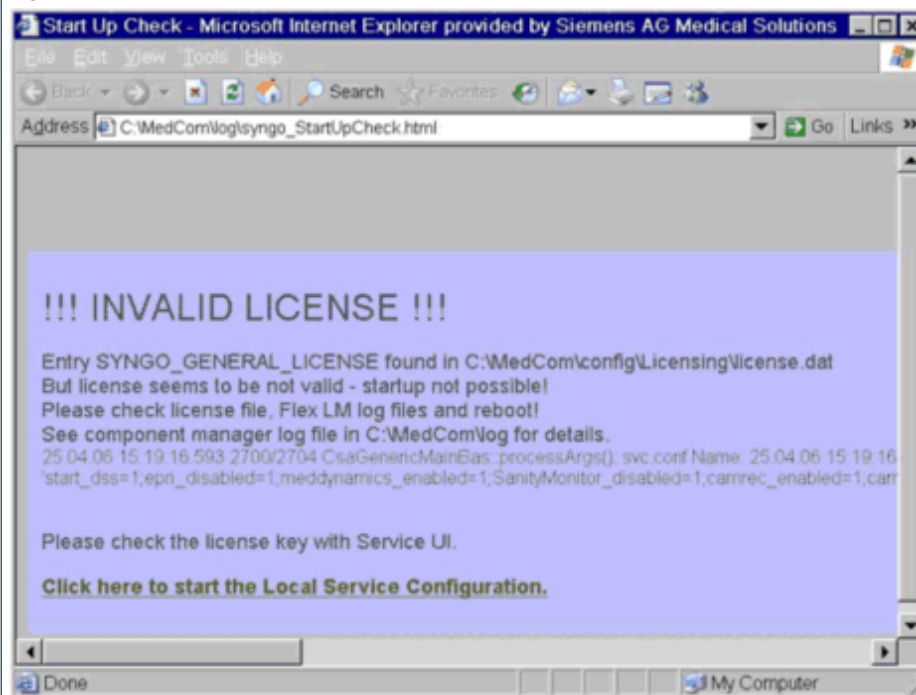

ACTION	MRC	MRSC
Start in local service by clicking the link "Click here to start the Local Service Configuration".	☐	☐
Use the new valid service key and click OK.	☐	☐
Select Configuration in the service platform.	☐	☐
Insert the floppy with the new license file	☐	☐
Select Licensing and import the new license file.	☐	☐
In the Licensing UI, choose "Select New License File" Fig. 120: 	☐	☐

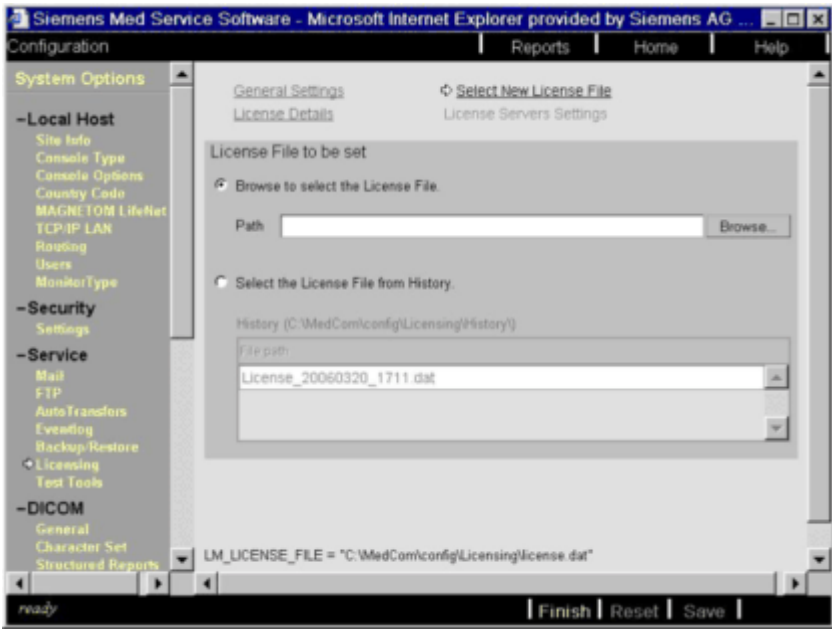ACTION	MRC	MRSC
- Under "Select New License File": Fig. 121:  - Use the Browse button to navigate to the source for the new license file - Import the new license file by clicking Save.	☐	☐
Go to Measurement "Meas. Settings" and press Save. Do not change any values on this page.	☐	n.a.
Reboot by clicking Home in the service UI.	☐	☐

5.1.5 Customer protocols

Only the **DotConfigurator** is supported for converting customer protocols within this software version.

5.1.6 Language settings and font size

Supported languages for meduser

- English (US)
- German (Germany)
- French (France)
- Spanish (Spain)
- Japanese
- Chinese (PRC)

Language change:


- Login as Meduser
- Activate "advanced user" (press Ctrl + Esc, type advanced User, select the program)
- Select: Options -> Configuration -> regional settings
Select the correct language in the format tab.
- Press OK for restart the system
- When the Syngo counting picture appears after reboot, activate the "advanced user" as soon as possible!
Wait until the syngo software is up again and check the language change.
- Select: Control Panel -> Appearance and Personalization.
Check if "Siemens_True_Color" is the selected theme.
- Run regedit (open registry) and check following key
HKLM\System\CurrentControlSet\services\eventlog\Application\<MessagedllPostfix e.g. AX_DR> "EventMessageFile" must contain a Messagedll with the correct language postfix (e.g. xxx_de.dll for german)
- Check the scanner symbol. if the MARS isn't connected the scanner symbol have to be crossed. Select this symbol and a new window with messages should appear


If you miss the point for activating the "advanced user" and the language change fail, don't panic. Switch back the language to the current format of the registry key and repeat the procedure above.



5.1.7 Keyboard and language setting

- Setting the keyboard layout and language

ACTION	MRC	MRSC
■ In the main UI, select Options -> Configuration.	☐	☐

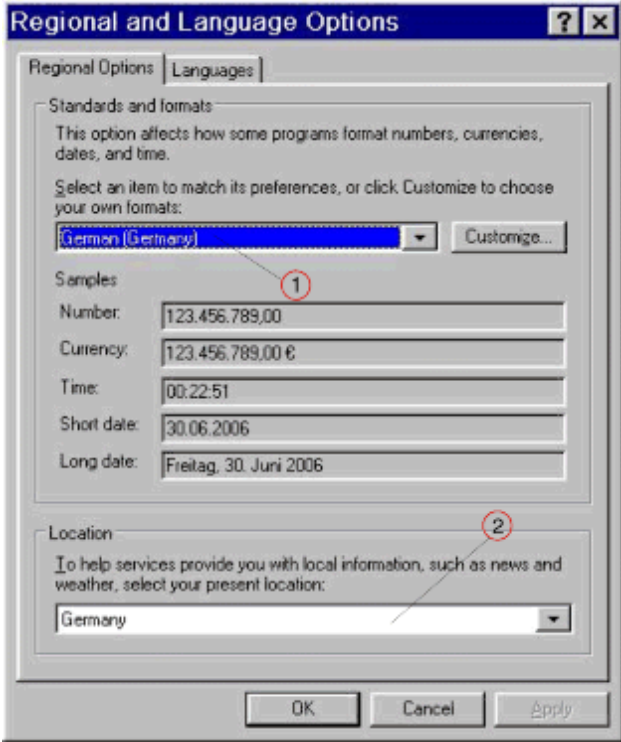ACTION	MRC	MRSC																
■ Click OK. Select the language for the main UI from the pull-down menu. The table below shows the supported entries for the languages. You can choose your country under "location". Fig. 122:  (1) Language (2) Country Only the following entries can be used for 1 in the graphic above: <table><tr><th>main UI language</th><th>entry in the pull down menu</th></tr><tr><td>german</td><td>German (Germany)</td></tr><tr><td>english</td><td>English (US)</td></tr><tr><td>english (metric)</td><td>English (UK)</td></tr><tr><td>french</td><td>French (France)</td></tr><tr><td>spanish</td><td>Spanish (Spain)</td></tr><tr><td>chinese</td><td>Chinese (PRC)</td></tr><tr><td>japanese</td><td>Japanese</td></tr></table>	main UI language	entry in the pull down menu	german	German (Germany)	english	English (US)	english (metric)	English (UK)	french	French (France)	spanish	Spanish (Spain)	chinese	Chinese (PRC)	japanese	Japanese	☐	☐
main UI language	entry in the pull down menu																	
german	German (Germany)																	
english	English (US)																	
english (metric)	English (UK)																	
french	French (France)																	
spanish	Spanish (Spain)																	
chinese	Chinese (PRC)																	
japanese	Japanese																	
■ Click Apply.	☐	☐																

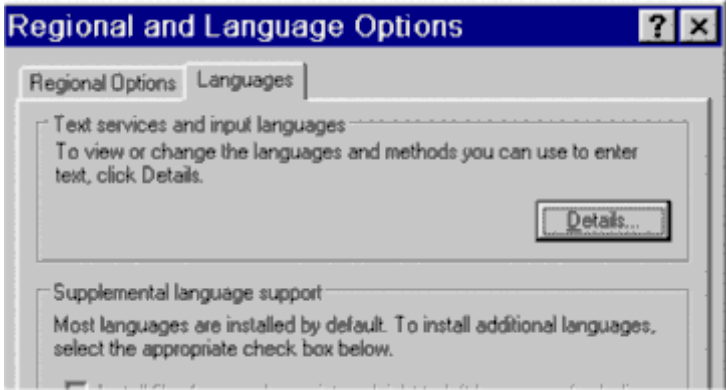
ACTION	MRC	MRSC
■ Select the Language tab in the Regional and Language Options task card. Fig. 123: 	☐	☐
■ Click Details.	☐	☐
■ In the "Text Services and Input Languages" window, press Add. Fig. 124: 	☐	☐

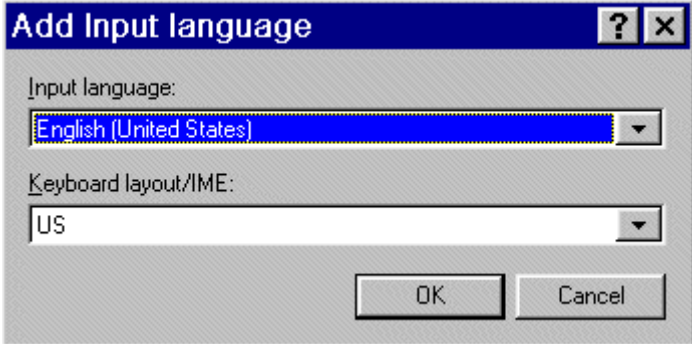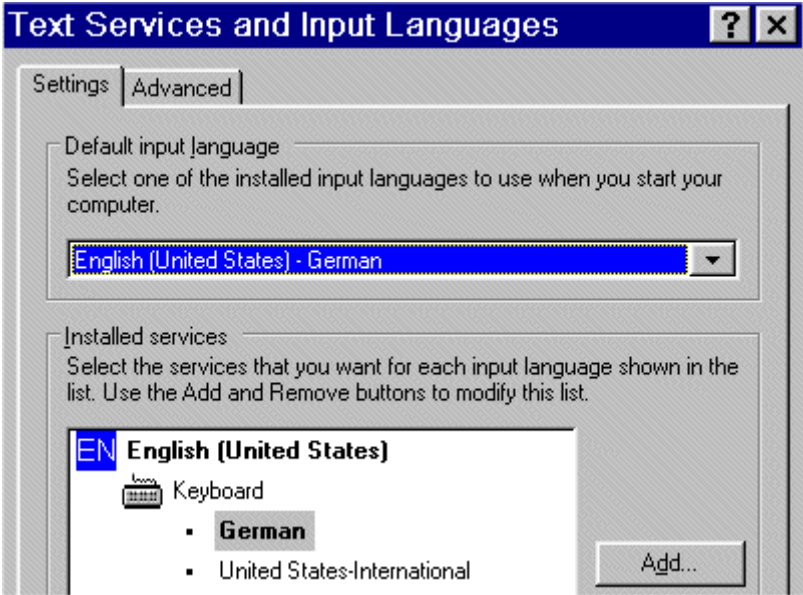ACTION	MRC	MRSC
- In the next pop-up window, choose the keyboard layout used under "Keyboard Layout." Select the language used for the main UI under "Input language". The input language also determines the system weight unit to be used. Fig. 125:  - English (US) uses pounds as the unit of weight and inches as the unit of length. - English (UK) uses kg as the unit of weight and cm as the unit of length.	☐	☐
Click OK.	☐	☐
Click Apply.	☐	☐
In the "Text Services and Input Languages" window, choose the newly configured keyboard layout under "default input language." Fig. 126: 	☐	☐
Click Apply	☐	☐

5.1.7.1 Keyboard and Language Settings - Japanese


- In the Service UI, select Options -> Configuration
- Select the DICOM Character Set > in the drop-down menu DICOM default character set, select **Japanese, Katakana**

Fig. 127: Configuration - DICOM default character set

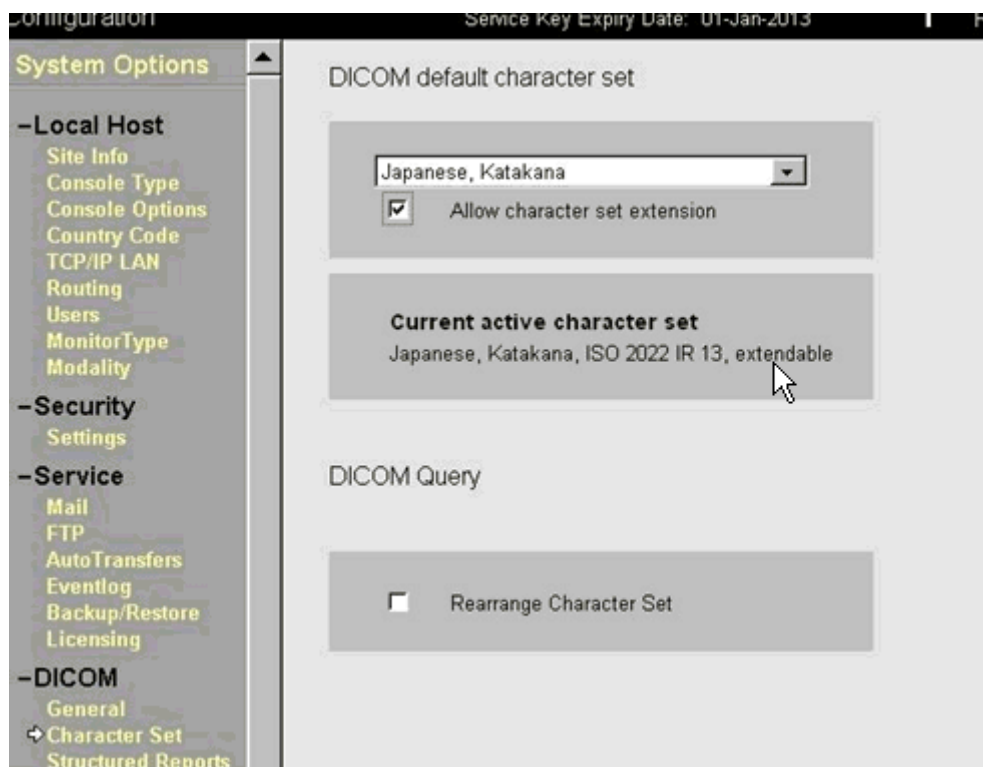


- Click **OK** to confirm the change in character set

Fig. 128: Configuration - character set data successfully saved

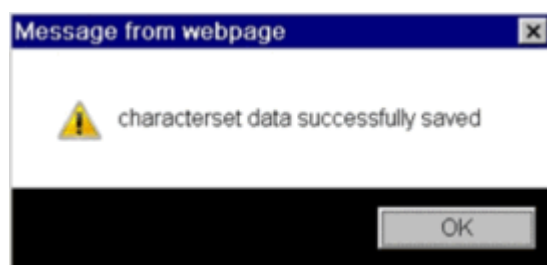


- Click **OK** to restart the application

Fig. 129: Configuration - a restart of the application SW is necessary to make changes valid

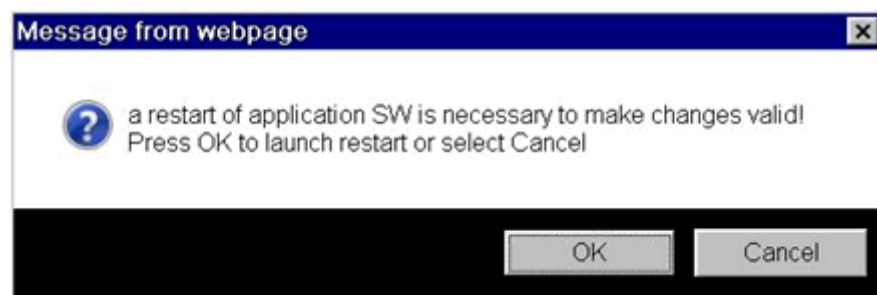


- Select **Options > Configuration**
- Select **Regional and Language Options**.

Fig. 130: Configuration panel

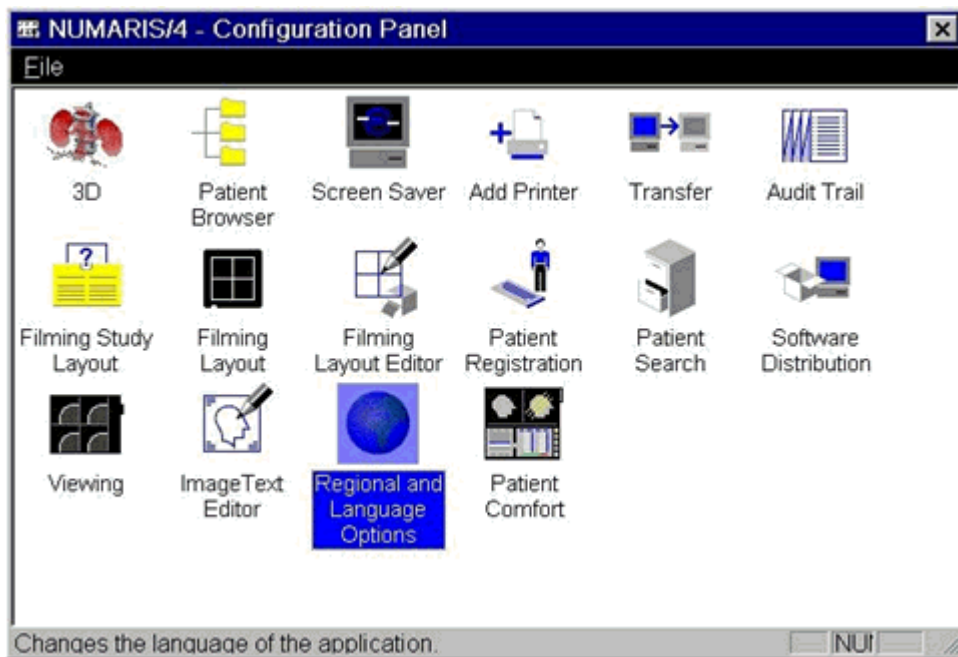

- Select **Japanese (Japan)** in the drop-down menu below “Format”
- Click **Additional settings**

Fig. 131: Regional and Language jp

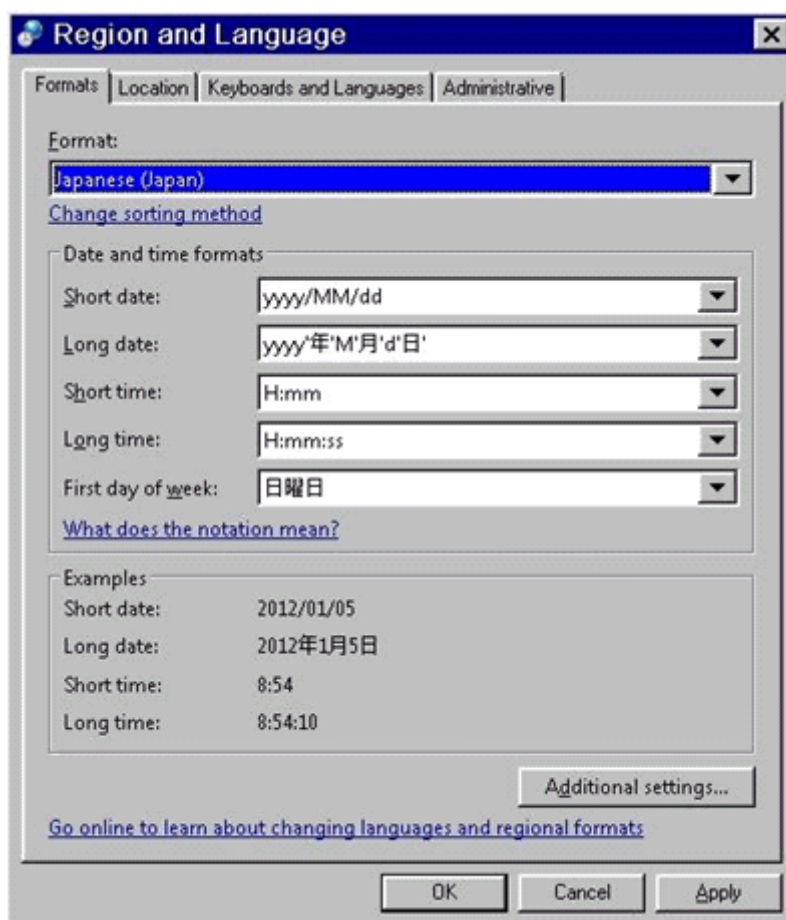


- Select the **Date** tab card
- Select the Calendar drop-down menu

Fig. 132: Calendar type Japan


和暦

Fig. 133: Customize Format - Calendar type

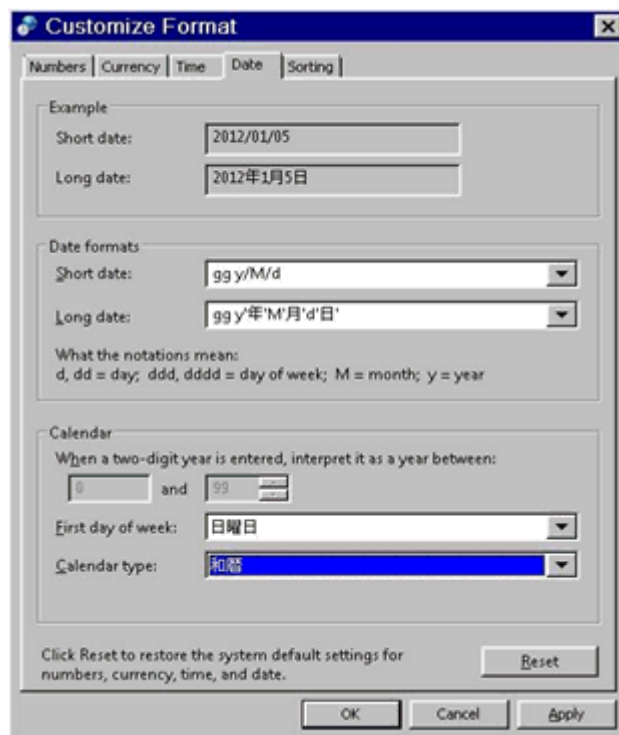

- In the drop-down menu "Current location:" select **Japan**
- Click **OK**

Fig. 134: Region and Language - Location Japan

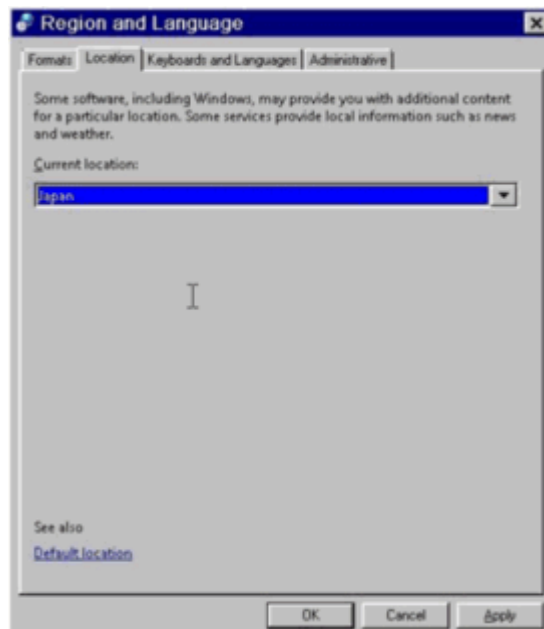


- Select the **Keyboards and Languages** tab card

- In the drop-down menu "Choose a display language:" select:

Fig. 135: Display language Japanese

日本語


Fig. 136: Region and Language - Keyboards and Languages Japan

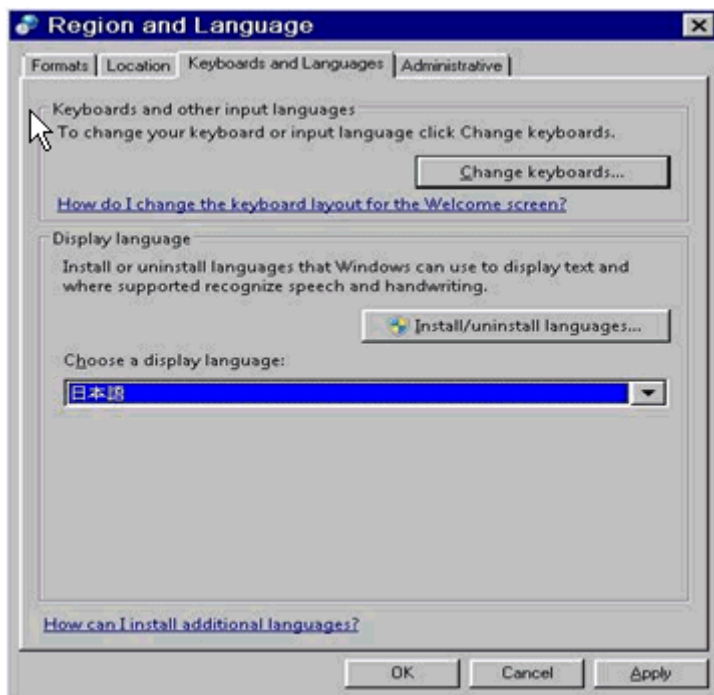

- Select the **General** tab card
- In the drop-down menu "Default input language:" select **Japanese (Japan) - Microsoft IME:**
- Add the Japanese keyboard if not available


- Select **Japanese (Japan) Keyboard Microsoft IME**

Fig. 137: Text Services and Input Languages - General Japan

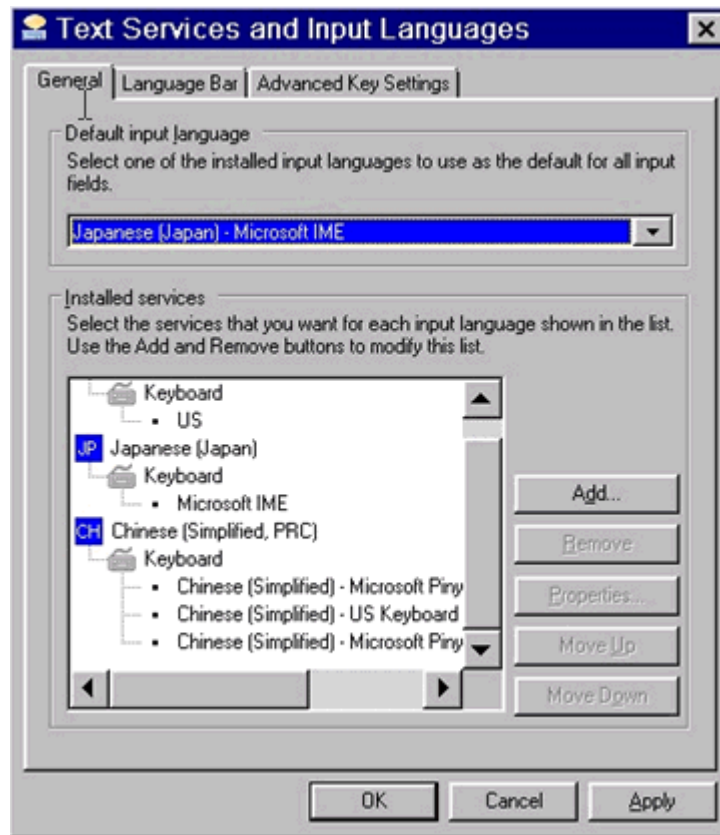

- Set the checkbox **IME** and click **OK**

Fig. 138: Add Input Language Japanese

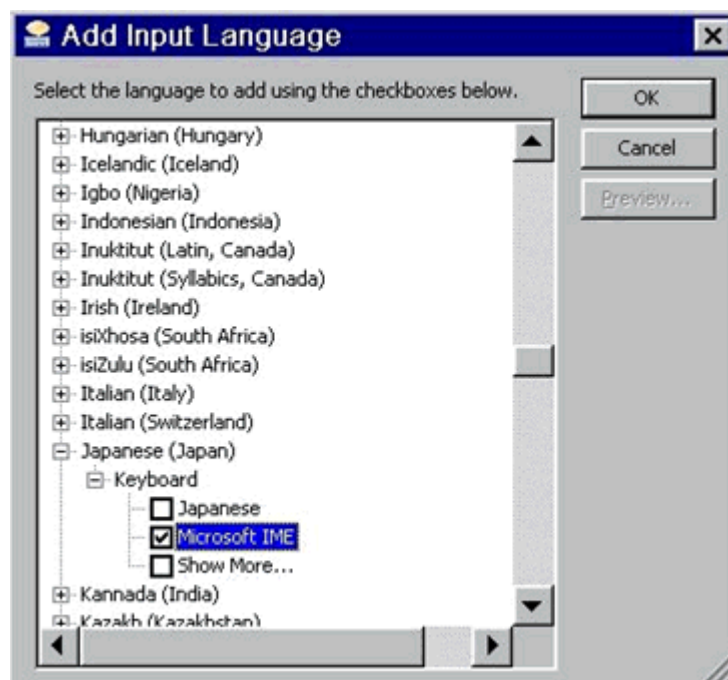

- Select the **Language Bar** tab card

- Select the **Floating On Desktop** radio button
- Set the checkboxes **Show additional Language bar icons in the taskbar** and **Show text labels on the Language bar**
- Click **OK**


Fig. 139: Text Services and Input Languages



- Select the **Administrative** tab card

- Click **Change system locale**

Fig. 140: Region and Languages - Administrative

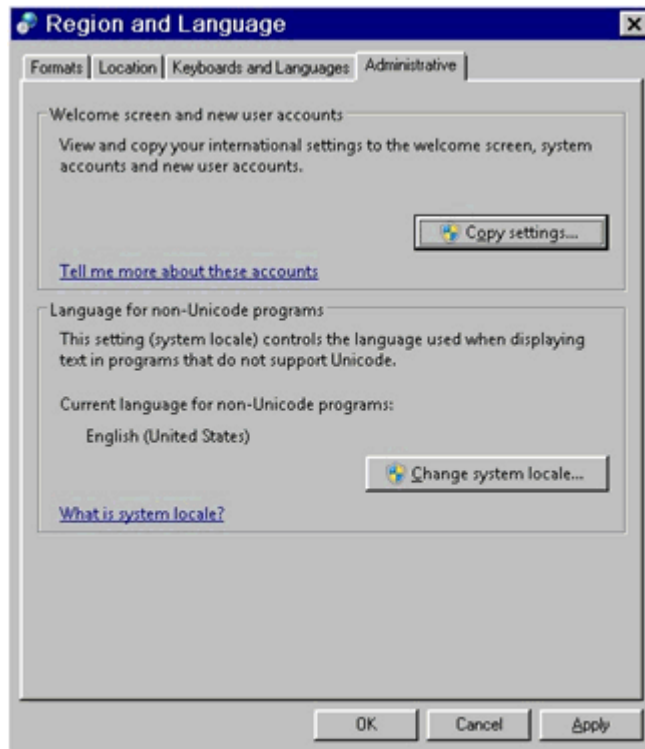

- Click **Apply** in the popup window

Fig. 141: Change Regional Options



- Select Japanese (Japan) in the next popup window


- Click **OK**

Fig. 142: Region and Language Settings - Japanese




- In the following popup window select **Restart now**

Fig. 143: Change System Locale - Restart now

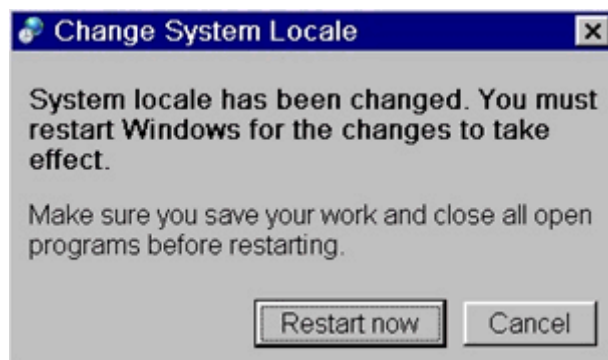


- Confirm this popup with **OK**

Fig. 144: OK - jp

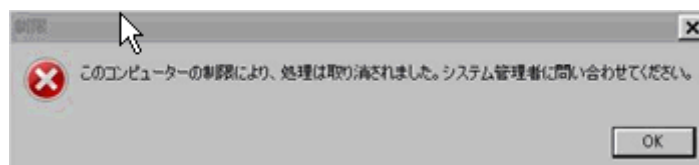

- Confirm the warning with **OK**


Fig. 145: Warning OK - jp

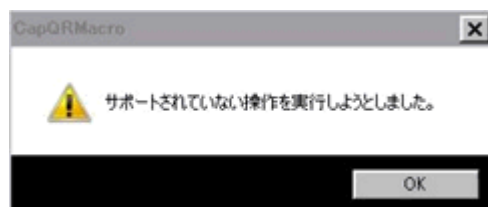

5.1.7.2 Switch to English from Japanese


- Select **Options > Service > Local Service**

Fig. 146: Options - Service - Local Service.jp

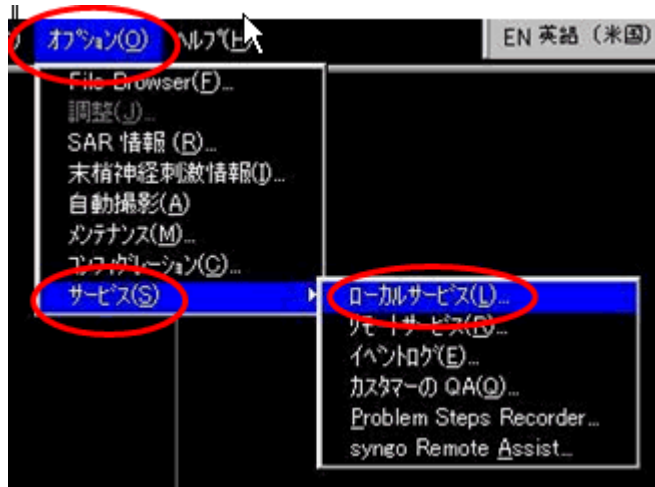

- Select **DICOM Character Set >** in the drop-down menu "DICOM default character set" select **Latin Alphabet No. 1**

Fig. 147: Configuration - Character Set - DICOM default character set (Latin Alphabet No 1).jp

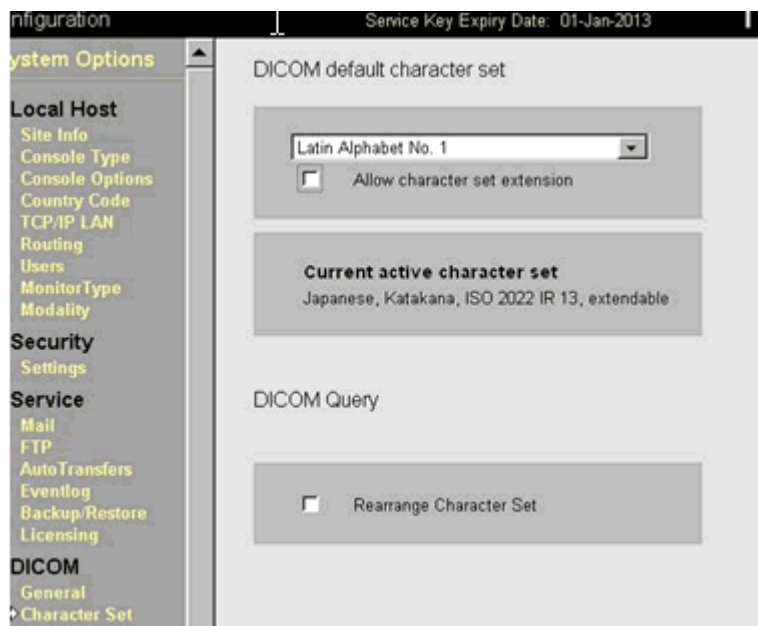


- Confirm “character set data successfully saved” popup with **OK**

Fig. 148: Character set data successfully saved - OK



- Confirm the “application restart popup” with **OK**

Fig. 149: A restart of the application SW is necessary - OK jp

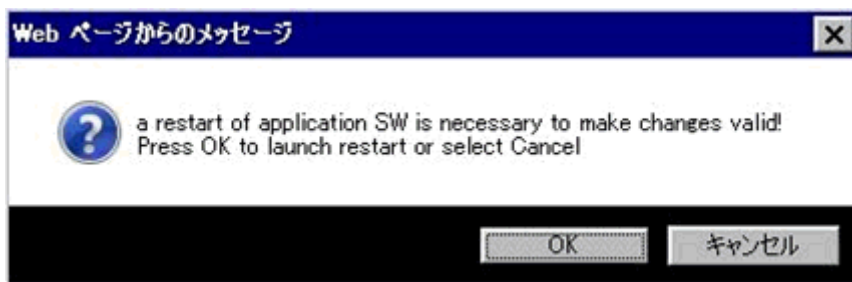


- When the application is up and running again, close the windows “Personalization” and “Local Service”

Fig. 150: Personalization and Local Service - jp




- Click **Options > Configuration**

Fig. 151: Options Configuration - jp

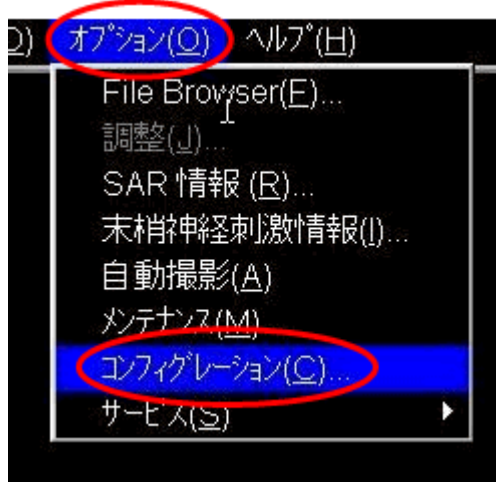

- In the Numaris/4 configuration panel, double-click **Regional and Language Options**


Fig. 152: Regional and Language Options - jp

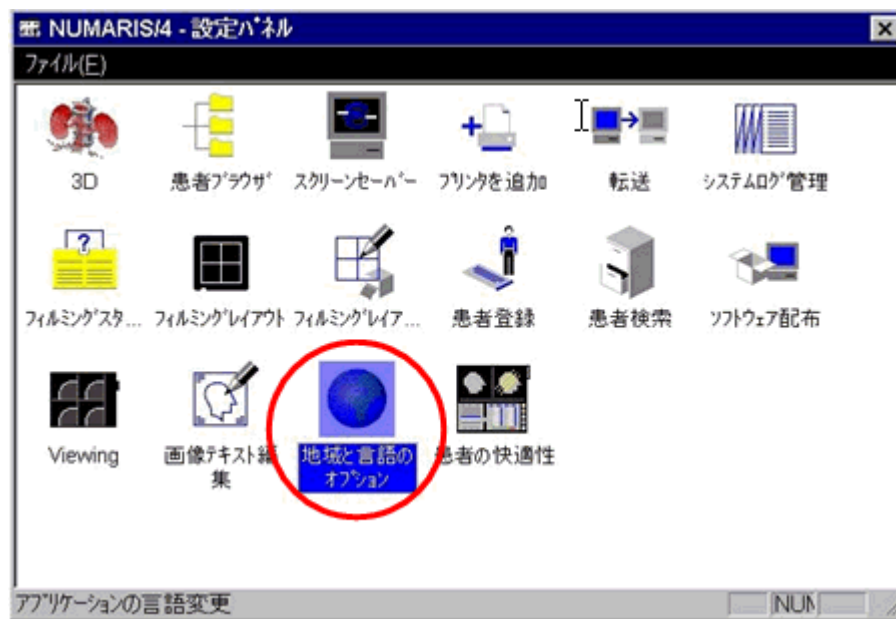

- Set the date time format as follows

Fig. 153: Region and Language - jp


英語 (米国)

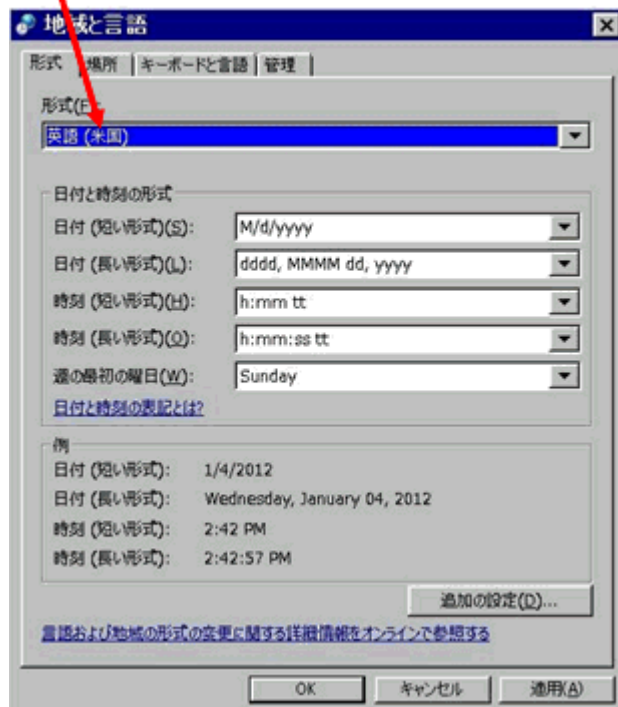

- Select the second tab


- Click **OK**

Fig. 154: Region and Language - Location jp

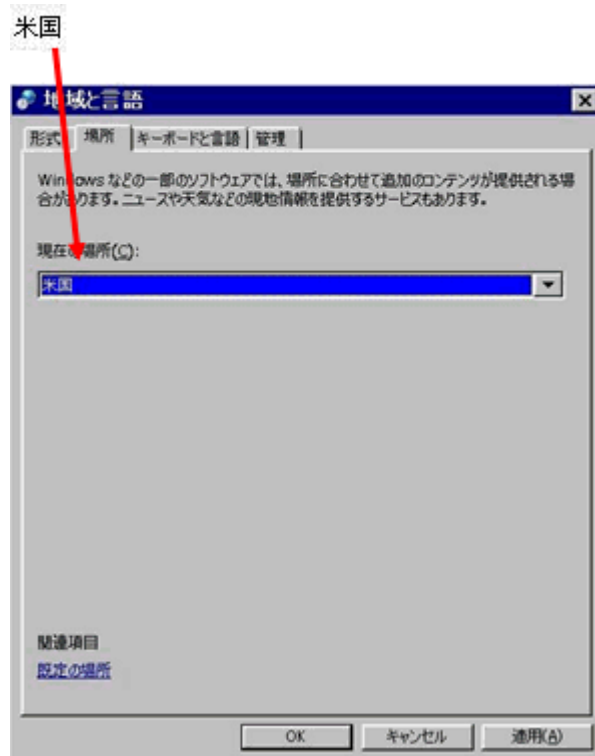

- Select **English**

Fig. 155: Region and Language - Keyboards and Languages jp

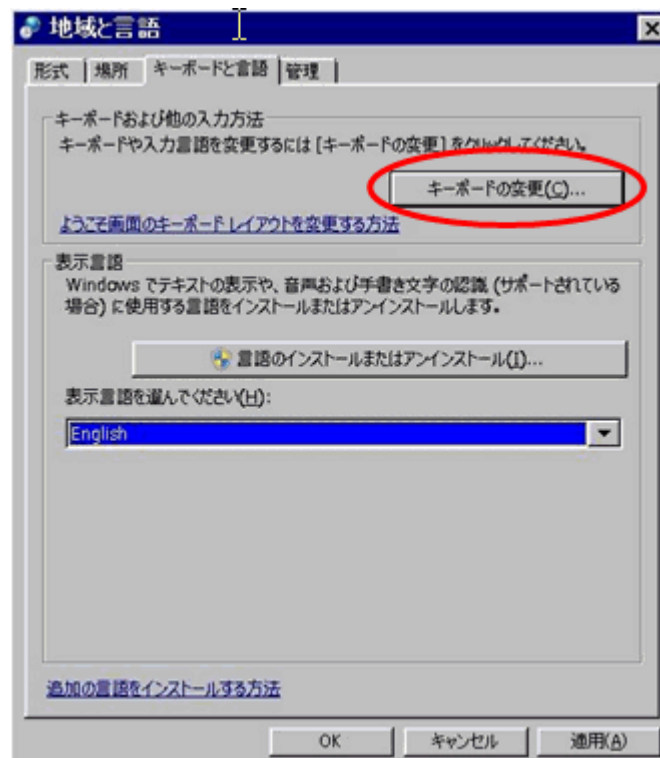


- Click the first tab from the left

- Select the following in the left drop-down menu


Fig. 156: Text Services and Input Languages jp

英語 (米国) - US

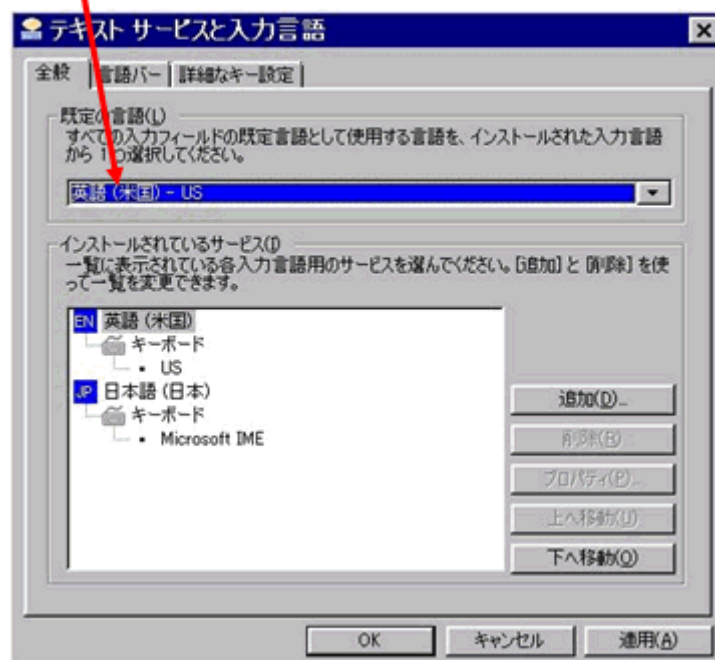

- Select the second tab > select the mid radio button
- Click OK

Fig. 157: Text Services and Input Languages jp

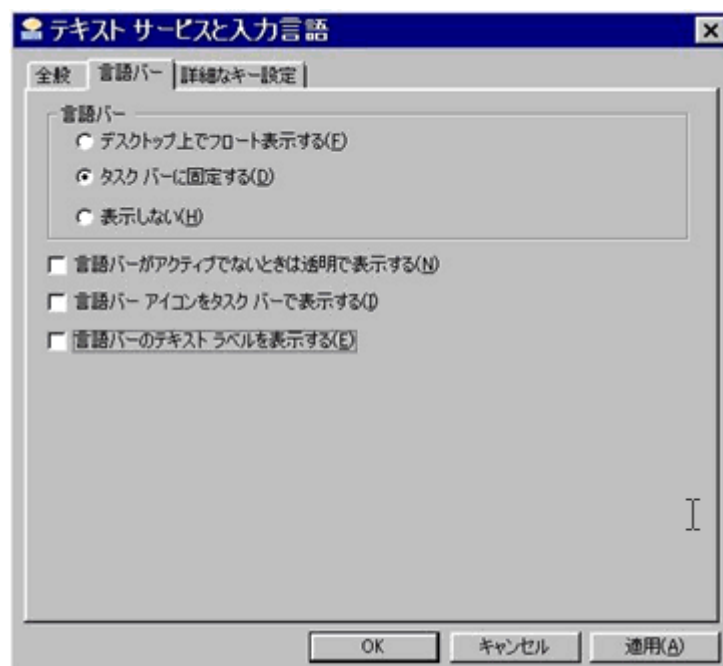

- Select the fourth tab > select the button

Fig. 158: Region and Language - Change system locale jp (zoomed)



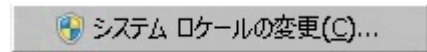


Fig. 159: Region and Language - Change system locale jp

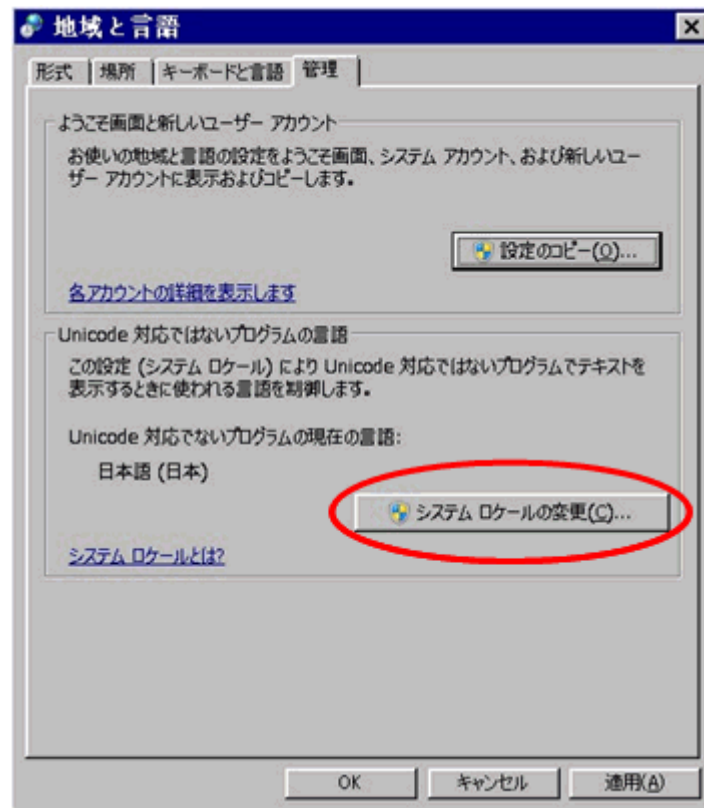


- Click Apply.

Fig. 160: Change Regional Options - Apply jp

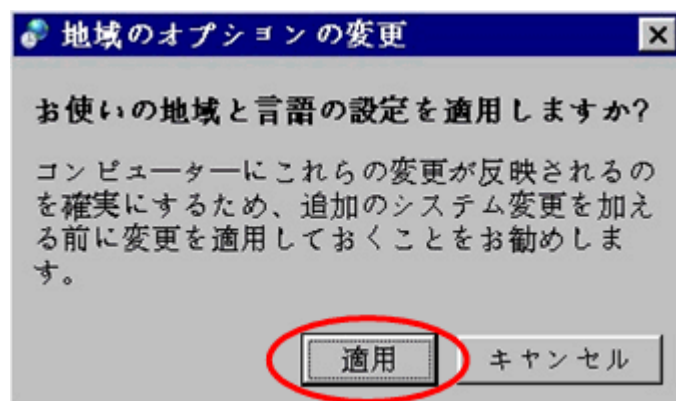


- Select the drop-down menu as follows and click **OK**

Fig. 161: Region and Language Settings - US English jp

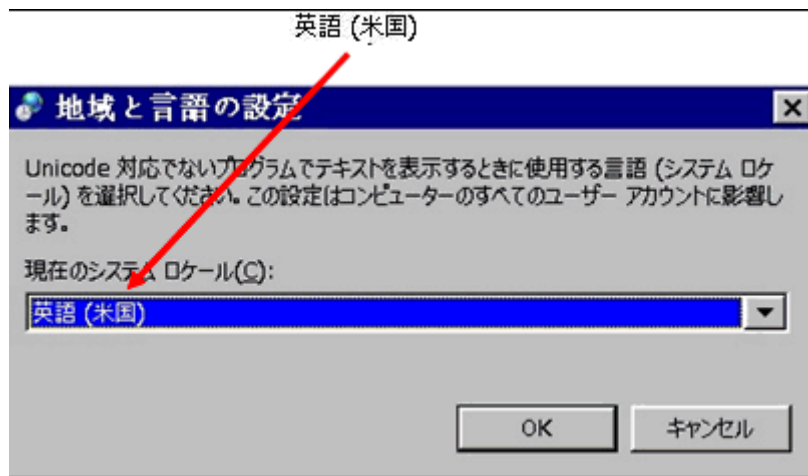

- Proceed as follows:

Fig. 162: Change System Locale - Restart now ch

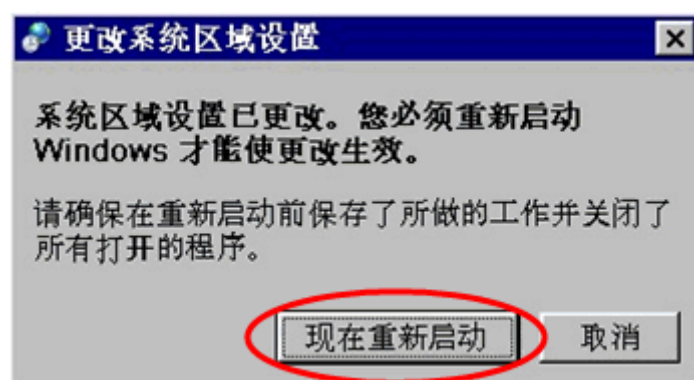

- Confirm the restriction popup with **OK**

Fig. 163: Confirm Restriction popup with OK



Changed from Japanese back to any other supported syngo UI language, the following occurs:


- After the reboot an error popup occurs. This error popup must be confirmed in order to start the syngo main UI.
- The screen settings configuration window is still opened and it must be closed manually.

5.1.7.3 Keyboard and Language Settings - Chinese

- In the Service UI, select **Options -> Configuration**


- Select the DICOM **Character Set** > in the drop-down menu "DICOM default character set" select **Chinese**

Fig. 164: Configuration - DICOM default character set

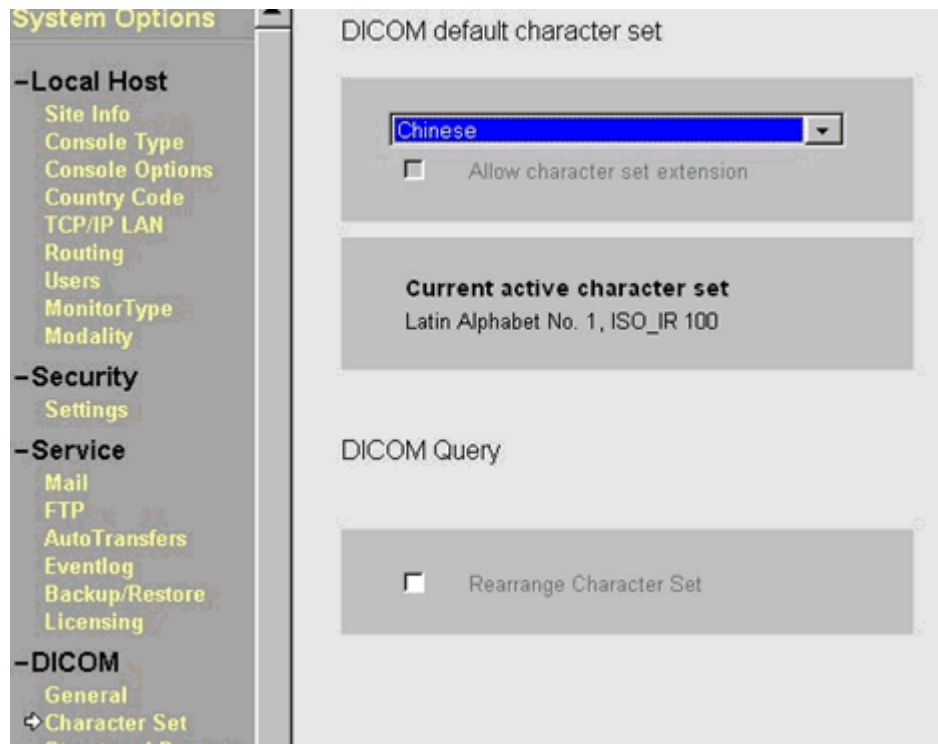


- Click **OK** to confirm the change in character set

Fig. 165: Configuration - character set data successfully saved

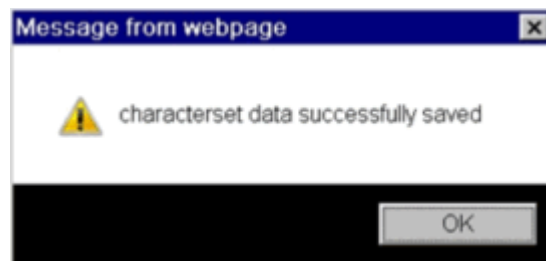

- Click **OK** to restart the application


Fig. 166: Configuration - a restart of the application SW is necessary to make changes valid

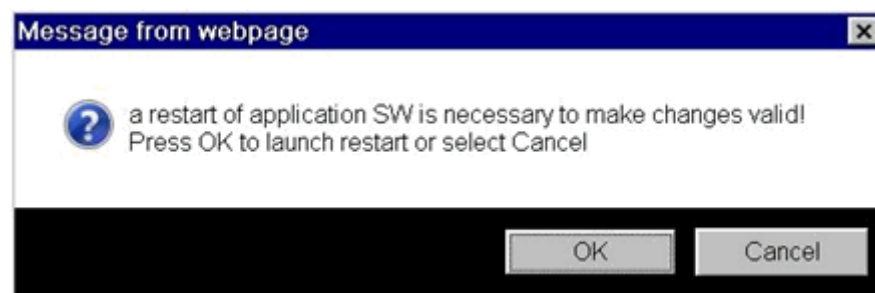

- Select **Options > Configuration**


- Select **Regional and Language Options**.

Fig. 167: Configuration panel

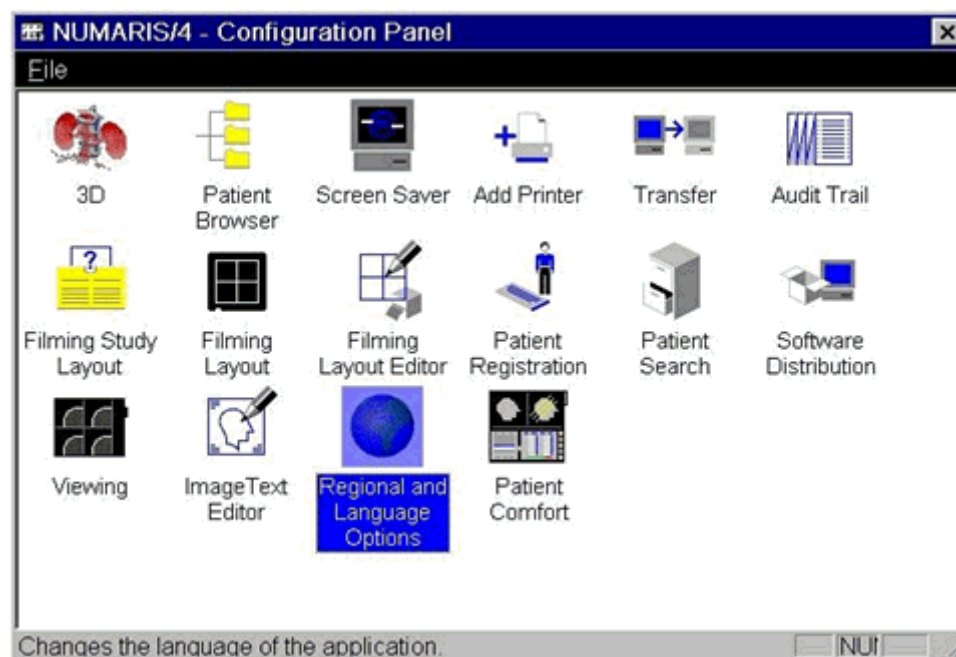

- Select **Chinese (Simplified, PRC)** in the drop-down menu below "Format"


Fig. 168: Regional and Language

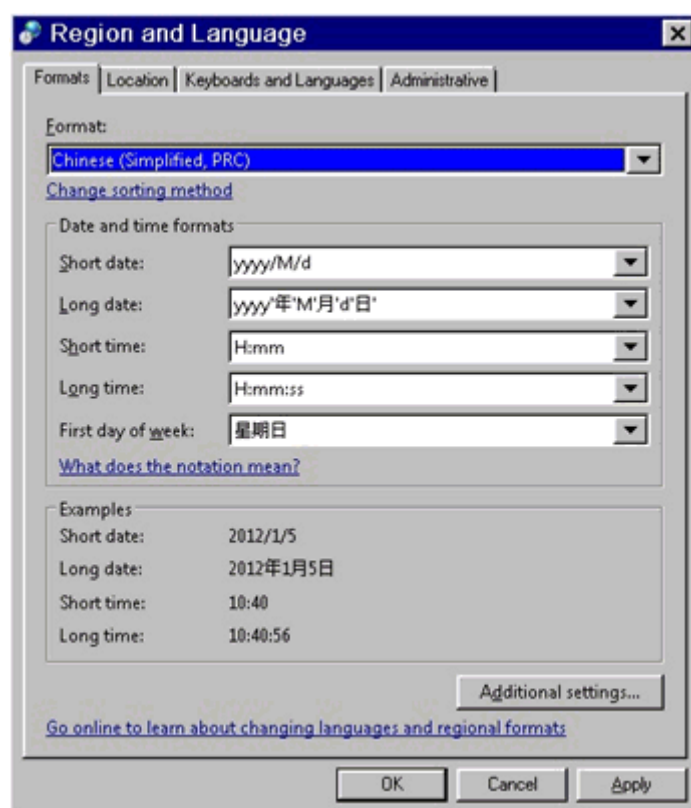

- Select the **Location** tab

- Select **Chinese** in the drop-down menu

Fig. 169: Region and Language - Location China

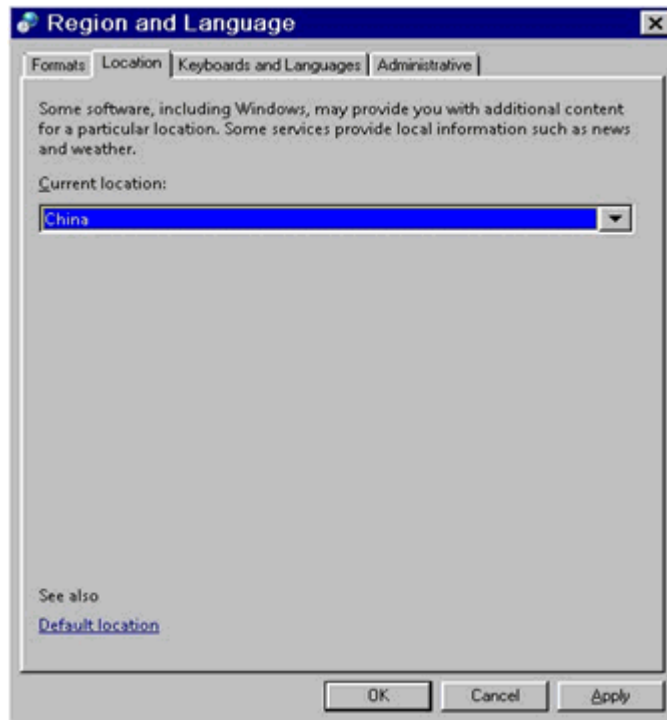


- Select the **Keyboards and Languages** tab
- Select

Fig. 170: Choose display language - Chinese

中文(简体)

- in the drop-down menu **Choose a display language** tab


Fig. 171: Region and Language - Keyboards and Languages

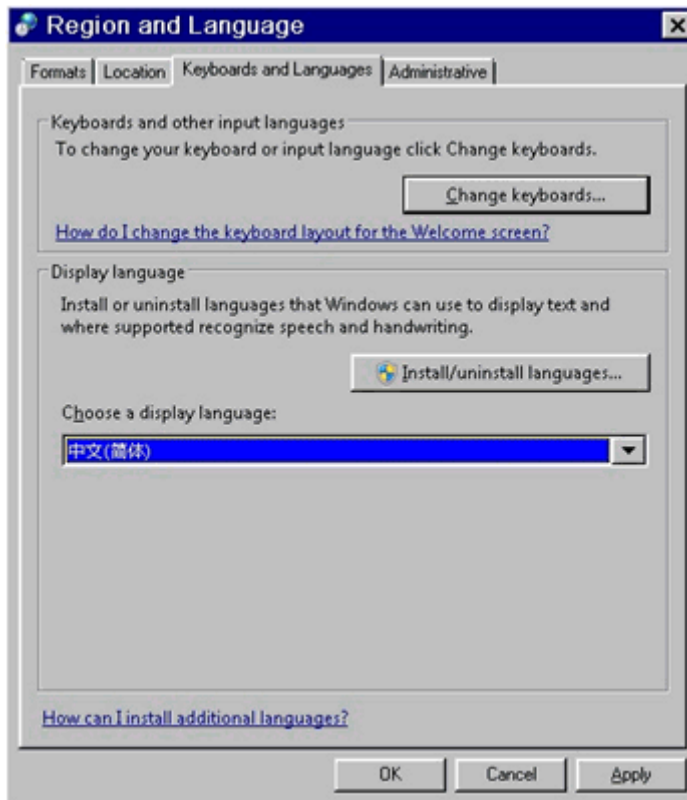

- Select the **General** tab.


- Select the drop-down menu **Chinese (Simplified, PRC) - Chinese (Simplified) - Microsoft Pinyin ABC Input Style**

Fig. 172: Text Services and Input Languages - General

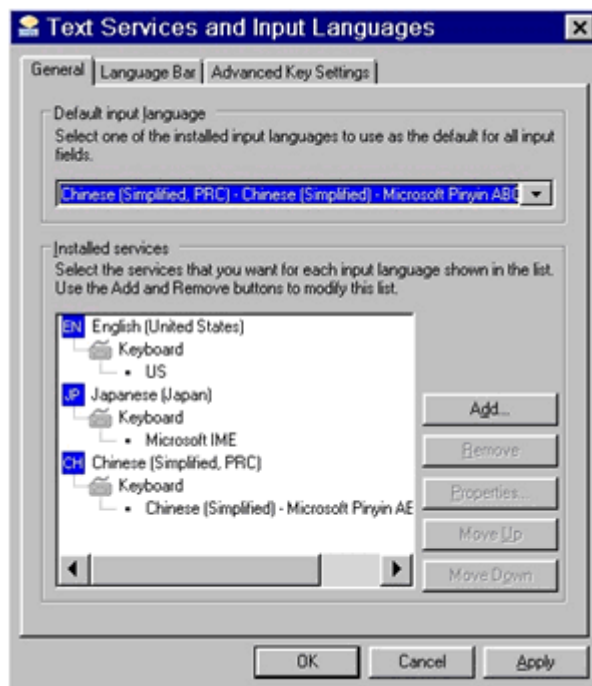

- Set the checkbox "Chinese (Simplified, PRC) - Chinese (Simplified) - Microsoft Pinyin ABC Input Style"


Fig. 173: Add Input Language

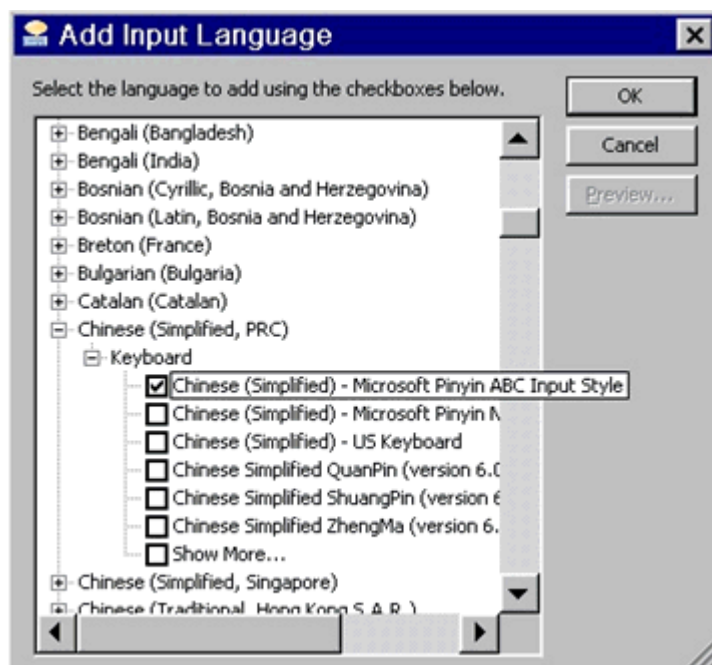

- Select the **Language Bar** tab card
- Select the **Floating On Desktop** radio button
- Set the checkbox **Show additional language bar icons in the task bar**

- Set the checkbox **Show text labels on the language bar**


Fig. 174: Text Services and Input Languages

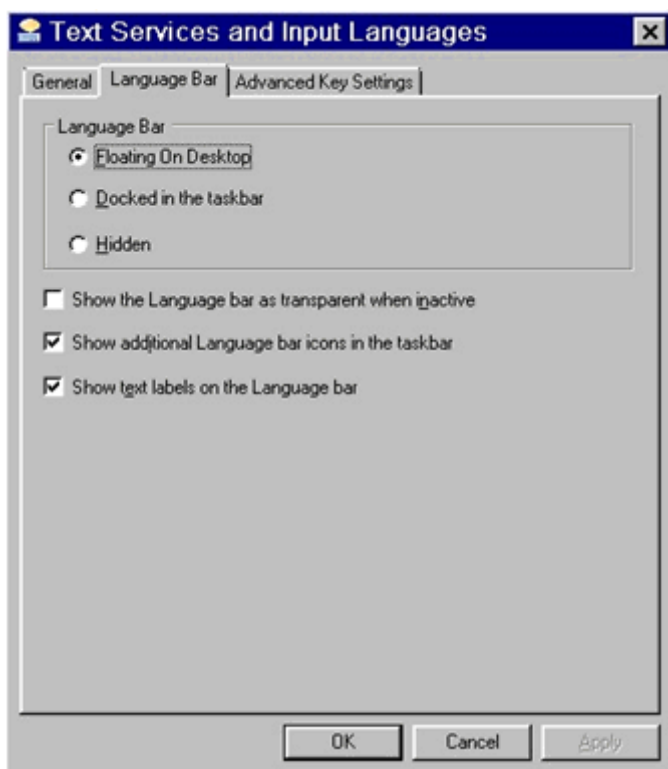

- Select the **Administrative** tab

- Click **Change system locale**

Fig. 175: Region and Languages - Administrative

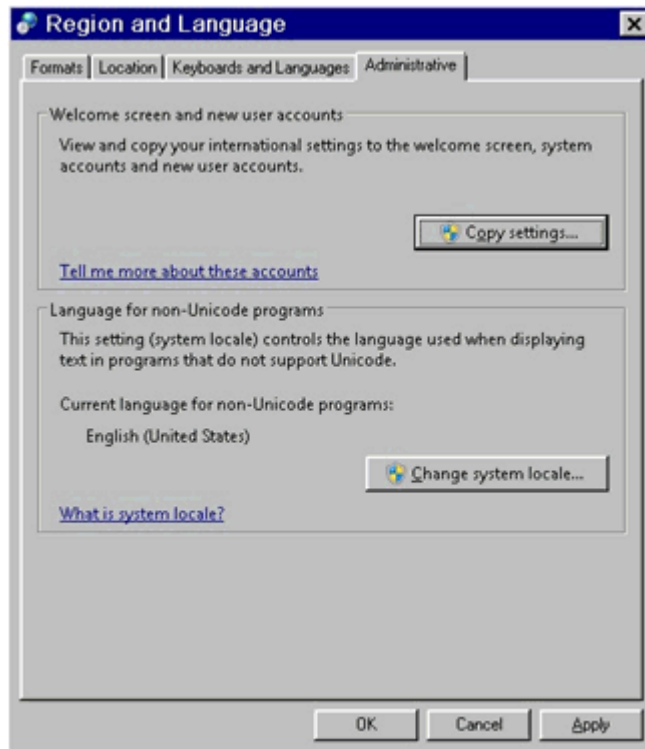

- Click **Apply** to confirm the region and language changes

Fig. 176: Change Regional Options




- Select **Chinese (Simplified, PRC)** in the drop-down menu “Current system locale”

Fig. 177: Region and Language Settings - Chinese




- Click **Restart now**

Fig. 178: Change System Locale - Restart now

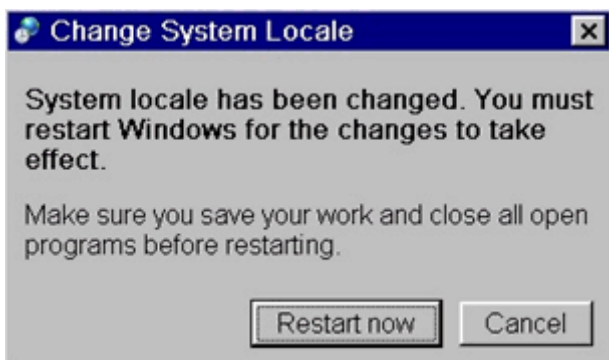

- Log in as meduser, wait until the application is up and running
- Confirm the following popups with **OK**

Fig. 179: OK - ch

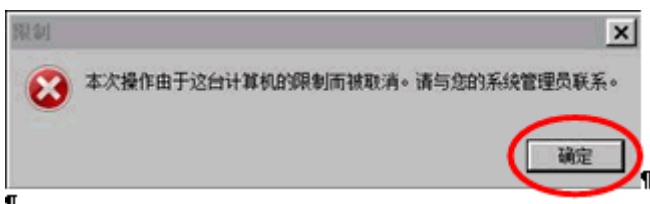

- Confirm the warning with **OK**

Fig. 180: Warning OK - ch




5.1.7.4 Switch to English from Chinese

- In the Service UI, select **Options -> Configuration**


- Select the DICOM **Character Set** > in the drop-down menu “DICOM default character set” select **Chinese**
- Open the Local Service service GUI

Fig. 181: Options - Service - Local Service ch

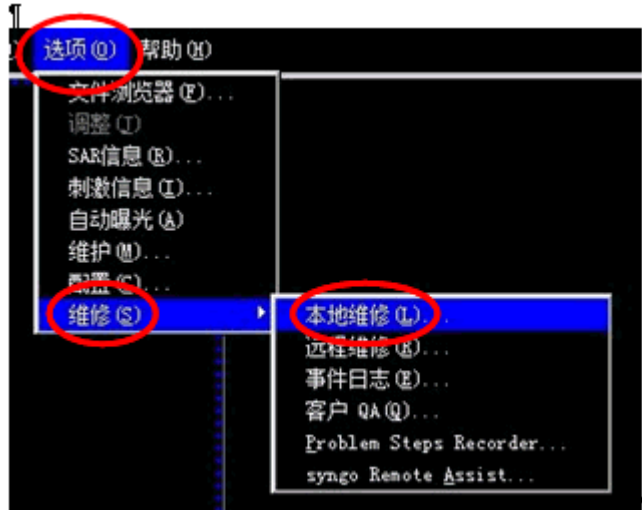


- Set the DICOM Character Set to **Latin Alphabet No. 1**

Fig. 182: Configuration - Character Set - DICOM default character set (Latin Alphabet No 1)

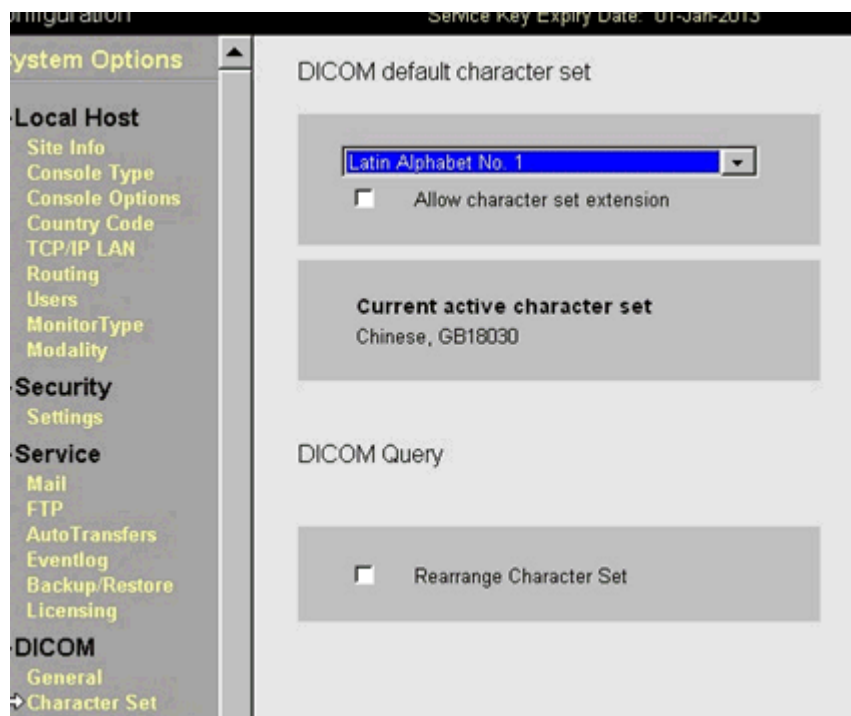

- Confirm the “character set data successfully saved” popup with **OK**

Fig. 183: Character set data successfully saved - OK ch

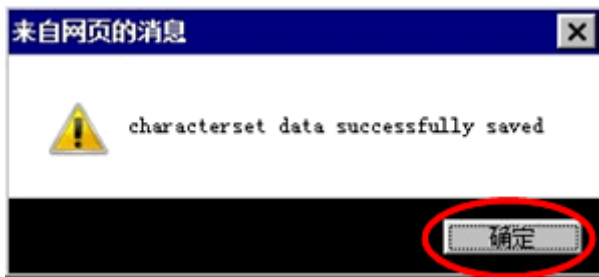

- Confirm the application restart popup with **OK**

Fig. 184: A restart of the application SW is necessary - OK ch




- When the application is up and running again, close the windows “Personalization” and “Local Service”

Fig. 185: Personalization and Local Service - ch



- Click **Options > Configuration**

Fig. 186: Options Configuration - ch

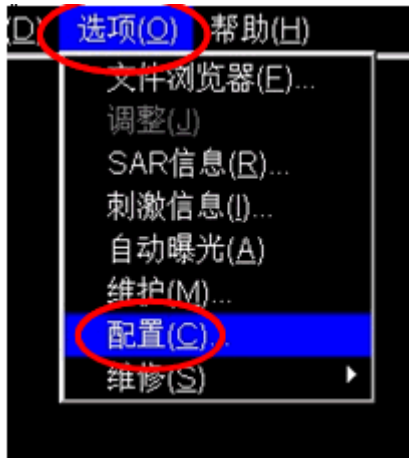


- In the Numaris/4 configuration panel, double-click **Regional and Language Options**

Fig. 187: Regional and Language Options - ch



- Set the date time format as follows


Fig. 188: Region and Language - ch

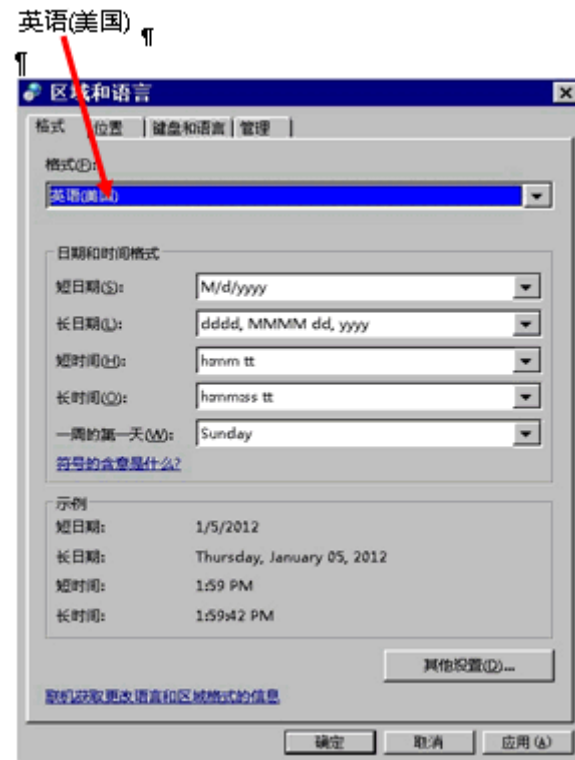

- Select the second tab


- Click OK.

Fig. 189: Region and Language - Location ch

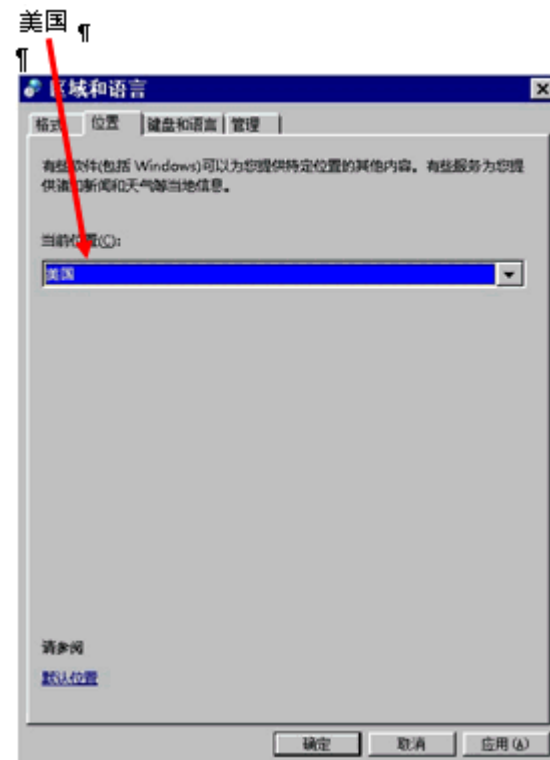

- Select **English**


Fig. 190: Region and Language - Keyboards and Languages ch

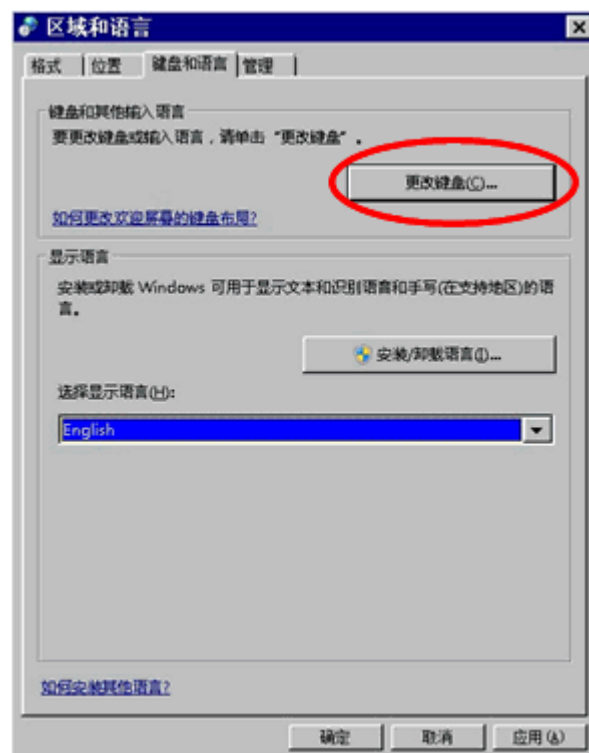

- Click the first tab from the left


- Select the following in the left drop-down menu

Fig. 191: Text Services and Input Languages ch

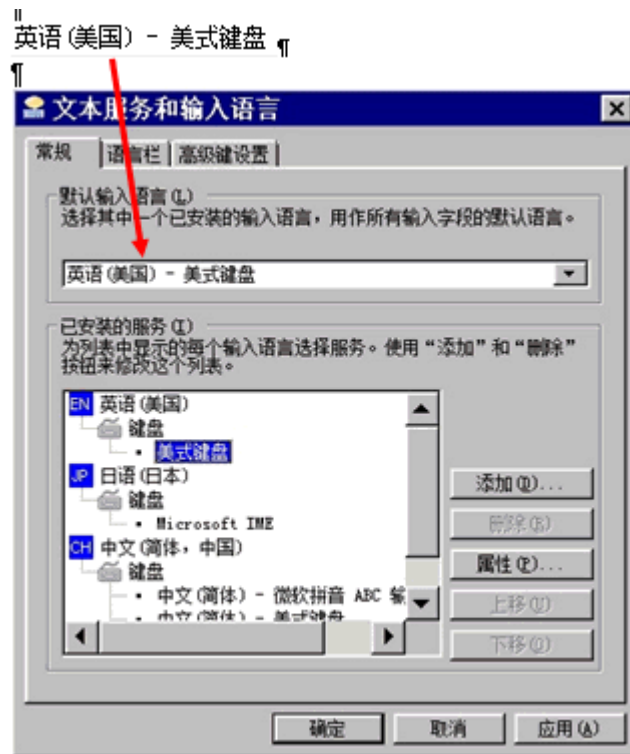

- Select the second tab > select the mid radio button
- Click OK

Fig. 192: Text Services and Input Languages ch

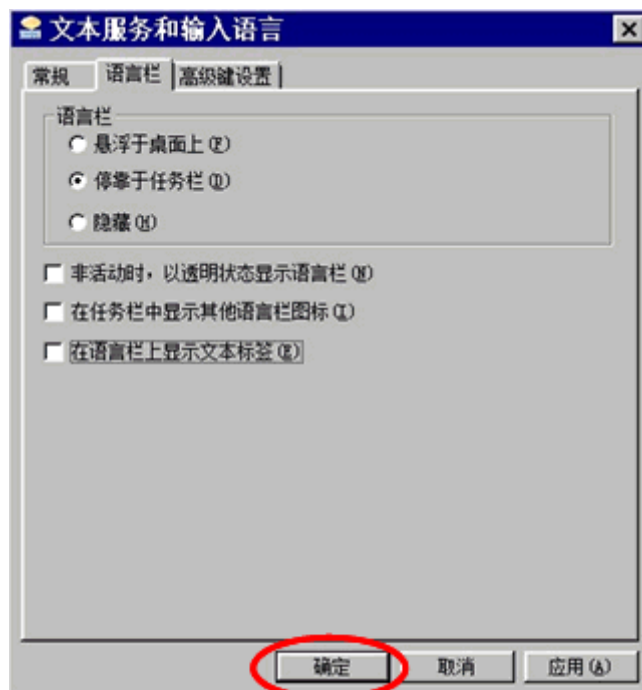

- Select the fourth tab > select the button

Fig. 193: Region and Language - Change system locale ch (zoomed)



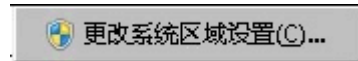


Fig. 194: Region and Language - Change system locale ch

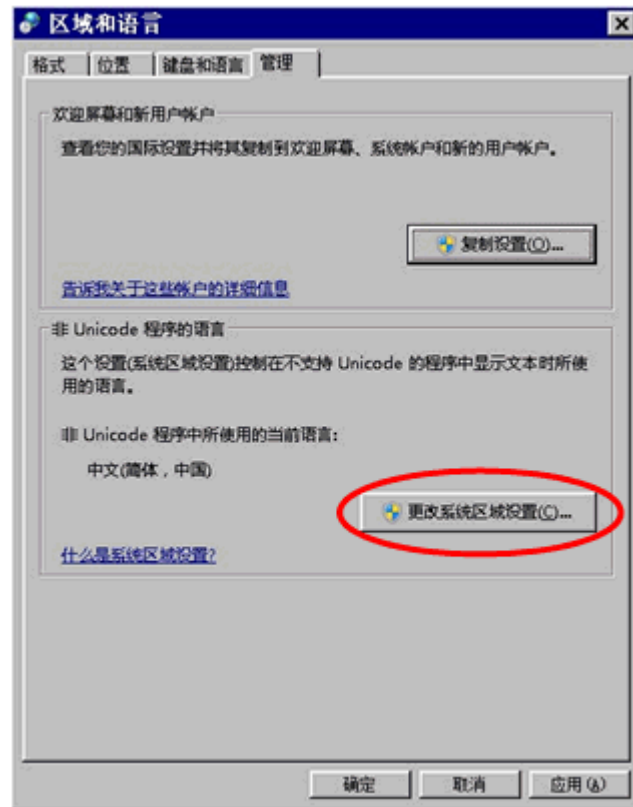


- Click Apply.

Fig. 195: Change Regional Options - Apply ch

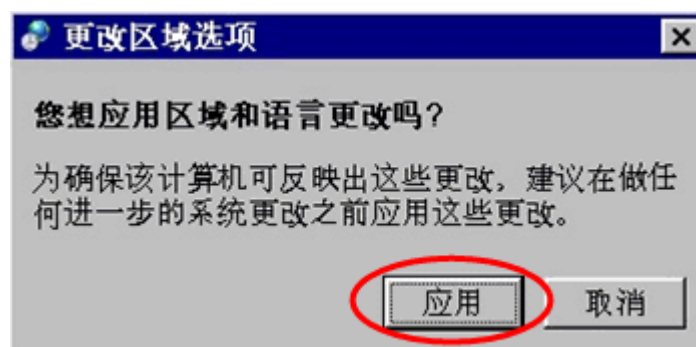


- Select the drop-down menu as follows and click **OK**

Fig. 196: Region and Language Settings - US English ch

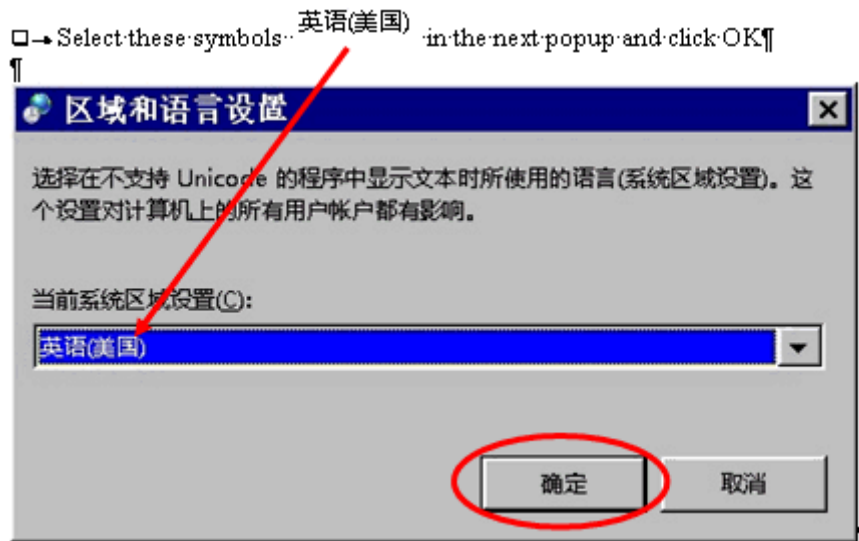

- Proceed as follows:

Fig. 197: Change System Locale - Restart now ch

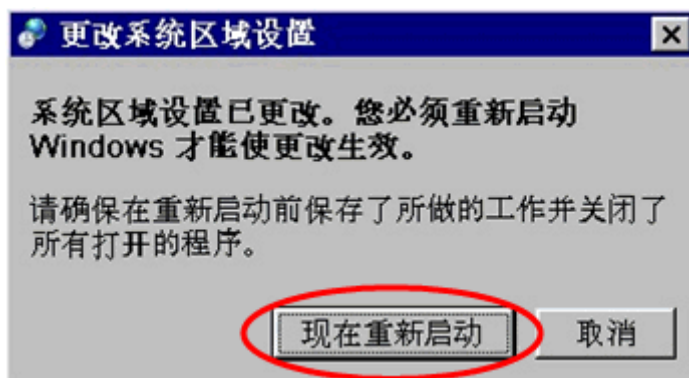

- Confirm the restriction popup with **OK**

Fig. 198: Confirm Restriction popup with OK



Changed from Chinese back to any other supported syngo UI language, the following occurs:


- After the reboot an error popup occurs. This error popup must be confirmed in order to start the syngo main UI.
- The screen settings configuration window is still opened and it must be closed manually.

5.1.8 Enable/Disable Inline Image Filter

The following procedure describes how to enable or disable the inline image filter

- Open local service.
- Select Configuration.
- Under Measurement select Meas. Settings
- In the Measurement Settings Page, go to the configuration entry 'Inline Image Filter' and select 'enabled' or 'disabled' in the pull-down menu.

Fig. 199: Service UI Measurement Settings

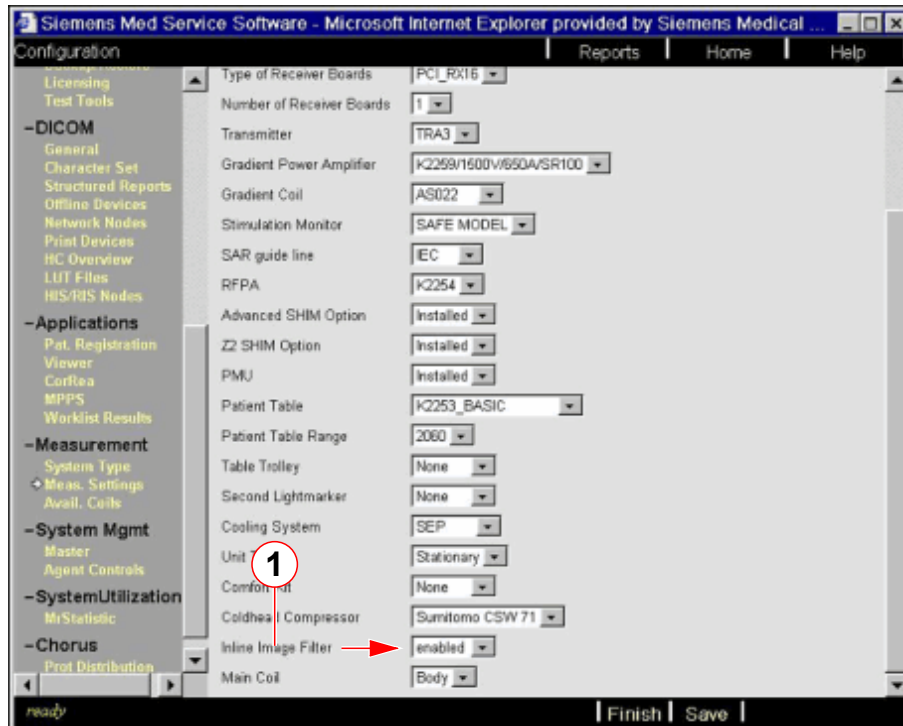

(1) Inline image filter

- After the changes have been made, press **Save** and wait until **ready** and **reboot system pending** are displayed in the status line of the service UI.
- Perform the reboot by pressing **Finish** and **Home**.
Confirm any displayed reboot pop-up windows.

5.1.9 Reset MRVCheck

The following error message starts to appear "Inform the local service, Finding occurred, check the log file c:\Medcom\logMrVCheck.log" during the host start-up process.

In this case, rename the file MrVcheck.ref under c:\medcom\config\mri\mrssystem to

MrVcheck.org

5.1.10 File cleanup for image creation

ACTION	MRAWP	MRWP
■ Login as administrator.	☐	☐
■ Start Windows Explorer: - Press <CTRL> + <ESC> to access the Start menu - In the Start menu, select All Programs ⇒ Accessories ⇒ Windows Explorer Windows Explorer starts.	☐	☐
■ Navigate to the directory C:\Medcom\temp\CDR_OFFLINE	☐	☐
■ Select all files and subdirectories in this directory and delete them.	☐	☐
■ Navigate to the directory C:\Medcom\MriDiagnostic	☐	☐
■ Select all files and subdirectories in this directory and delete them.	☐	☐
■ Navigate to the directory C:\Medcom\MrLogRepository	☐	☐
■ Select all files and subdirectories in this directory and delete them.	☐	☐
■ Navigate to the directory C:\Medcom\MrLogRepository	☐	☐
■ Select all files and subdirectories in this directory and delete them.	☐	☐
■ Navigate to the directory C:\Temp\HostLog	☐	☐
■ Select all files and subdirectories in this directory and delete them.	☐	☐
■ Navigate to the directory C:\Users\meduser\AppData\Local\Microsoft\Windows\WER\ReportQueue	☐	☐
■ Select all files and subdirectories in this directory and delete them.	☐	☐
■ Empty the Recycle bin: - Select the icon Recycle Bin on the desktop - press the right mouse button - in the displayed menu select the entry Empty Recycle Bin	☐	☐
■ Close the windows explorer	☐	☐

5.2 Hints

5.2.1 Exporting and backing up customer measurement programs

- Exporting customer measurement programs
 - **Open/Start** Exam Explorer.
 - In the user tree, select the region to be stored/saved.
 - Select the **Object** menu item, then select **Export**, and select the path C:\ExamDBBackup. Confirm by clicking the **Export** button.
 - Repeat steps 2 and 3 for all selected regions.
 - Close Exam Explorer.
- Making a Backup
 - Select the **Options** menu item, then select **Local Service**.
 - Enter the valid service key and confirm by clicking **OK**.
 - Select **Backup/Restore**.
 - Use the **Backup** command to select the CD writer and the **ExamDB** package.
 - Click **Go** to save the customer measurement programs to the CD-R drive.

5.2.2 Restoring and importing customer measurement programs

- Restoring the backup
 - After completing the installation of the MRease Software, select the **Options** menu item, then select **Local Service**.
 - Enter the valid service key and confirm by clicking **OK**.
 - Select **Backup/Restore**.
 - Use the **Restore** command to select the CD writer and the **ExamDB** package.
 - Click **GO** to restore the programs.
- Importing customer measurement programs
 - Open Exam Explorer.
 - Select the user tree.
 - Select the **Object** menu item, then **Import**, and select path C:\ExamDBBackup. Confirm the region to be imported by clicking the **Import** button.
 - Repeat step 3 for all selected regions.
 - Close Exam Explorer.

5.2.3 Information for the account administrator

Do not mark the checkbox "Random Password" under local service, configuration, Users section Account for Administrator. If this is marked, it is no longer possible to log in as the administrator because the password is no longer known.

5.2.4 Changing network addresses for the measurement network

If the customer uses the same network address range as the measurement network, the addresses for the measurement network can now be changed according to the following procedure:

- Start Local Service
 - Go to Configuration
 - Select TCP/IP
1. Select the network adapter
`[00000001] Measurement (internal) Network (network to the imager)`
 2. Select the IP address, click the **Del** button in order to remove the entry.
 3. Add the new IP address by inserting it into the appropriate field and subnet mask. Click the **Add** button.
 4. Click **SAVE** to save the changes.
 - Reboot the MRC
 - Wait until the main UI is displayed.
 - Shut down the system and switch it off
 - Switch on the system.

Now the system should come up with the new measurement network IP addresses.

5.2.5 Use of customer-built coils or 3rd party coils with Magnetom systems

5.2.5.1 Definition:

Customer coils: Non-Siemens product coils, coils not from the accessory list

5.2.5.2 General:

A small number of research-oriented customers has a special contract with Siemens about customer coils interfacing with the system. These customers received specific coil interface information and a site-specific password to connect customer coils at their system.

There are also customers who work with coils from the coil manufacturer RAPID. RAPID delivers customer coils with valid coil files with and without compatibility declaration.

5.2.5.3 Coil codes:

Siemens has reserved coil codes for customer coils: 93, c2, a2, a3, a1, 91, 71, 81

5.2.5.4 Customer coil file directory:

The customer coil files have to be stored in the directory C:\MedCom\MriCustomer\Coils

5.2.5.5 Template coil files in the Siemens coil file directory:

Template coil files for customer codes 93, c2, a2, a3 are stored in the Siemens coil directory C:\MedCom\MriSiteData\measurement\Coils. If one or more of these codes is used in the customer coil file directory, the respective template files have to be renamed in the Siemens coil file directory.

5.2.5.6 Tip for SW update:

In case of a SW update, the template coil files for customer codes 93, c2, a2 and a3 are installed again in the Siemens coil file directory. If one or more of these codes are used in the customer coil file directory, the respective template files have to be renamed again in the Siemens coil file directory.

5.2.6 Startup error "Found Databases with Equal IDs" on MRSC

5.2.6.1 Description:

After the restore of the dynamic data, the MRSC with shared database is not starting up. During startup, the following error message is displayed:

```
ERROR: Found Databases with Equal IDs: xx
startup_error_code=CSA_DB 0xE800004B
```

5.2.6.2 Solution:

Check the host name of the MRC and MRSC. The host name should not contain the string **ID** or **id**. If the host names on MRC/MRSC containing the string **ID** or **id**, the host names must be modified. Otherwise, the start-up check fails again.

5.2.7 Drive letter assignment for the CD/DVD ROM drive and/or for the CD/DVD-R drive

5.2.7.1 Description:

The drive letter assignment for the CD/DVD ROM drive and/or for the CD/DVD-R drive is not correct. The following is displayed in the local service configuration

Fig. 200: Error message



Fig. 201: Error message


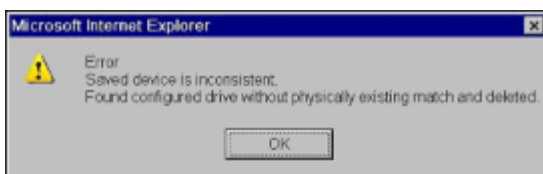


Fig. 202: Local service UI

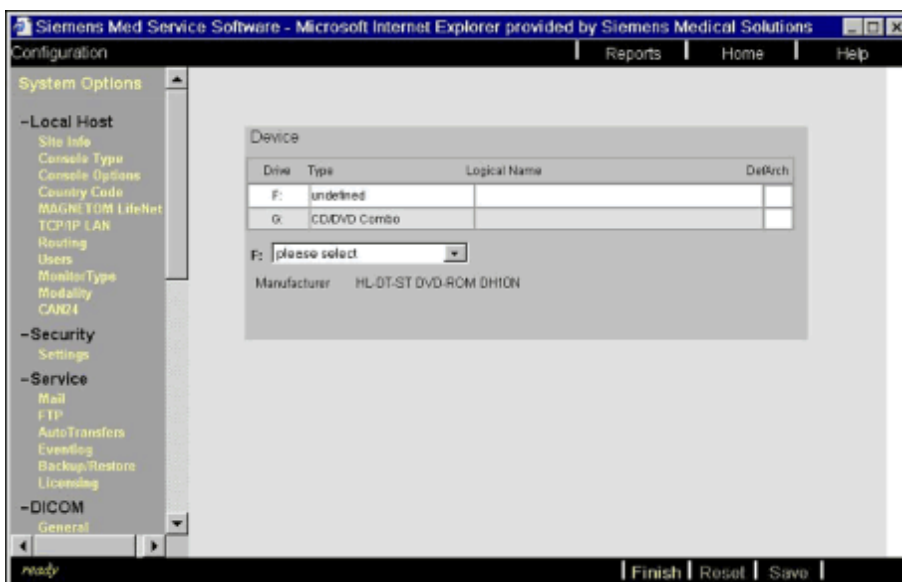


5.2.7.2 Solution:

The following procedure describes the correction for the drive letter assignment.

- Log in as administrator
- Start the disk manager by selecting **START --> RUN** and entering `diskmgmt.msc`


Fig. 203: RUN pop-up

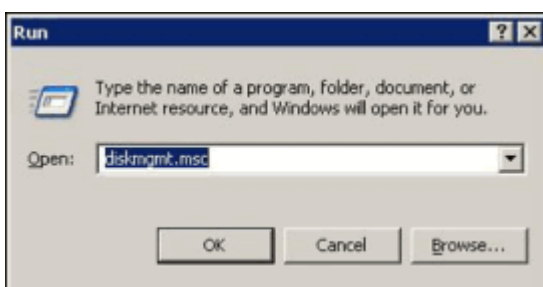

- The disk manager pops up.
- Select the drive to be (re-) assigned with a right mouse click, select '**Change drive letter and paths...**' from the pop-up menu, select '**Change**' in the next pop-up and choose the proper drive letter from the drop-down list of the next one. Click '**Yes**' in


the 'Confirm' pop-up. (Repeat this sequence for every drive for which you have to change the drive letter).

Fig. 204: Drive letters are wrong

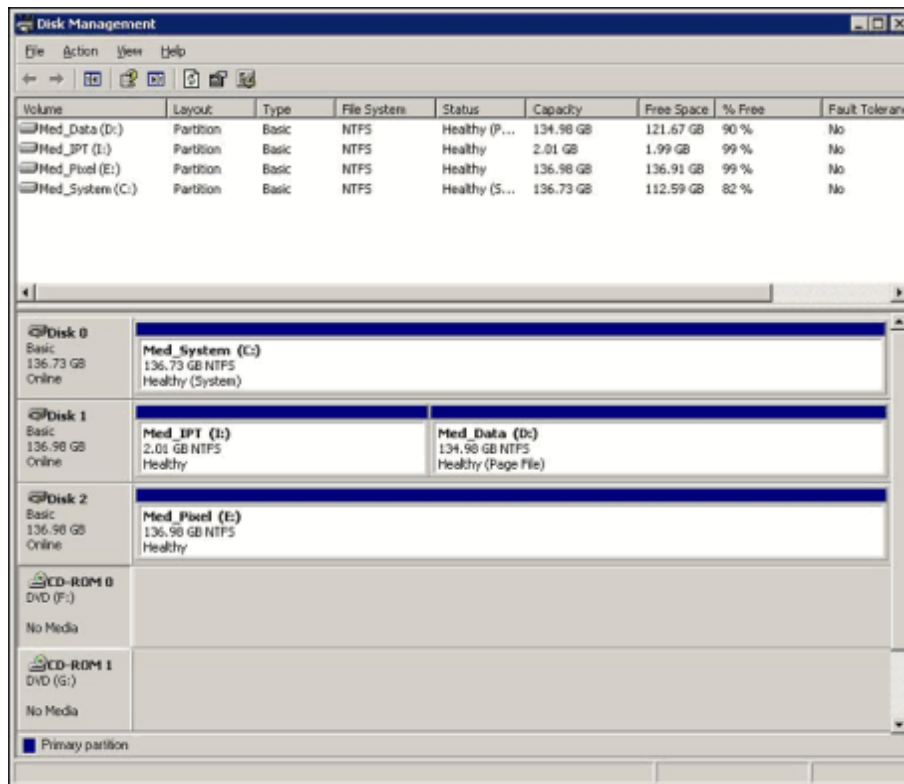

- CD/DVD-ROM drive -> **R:**
- CD/DVD-R (writer) drive-> **S:**

Fig. 205: Drive letters are OK

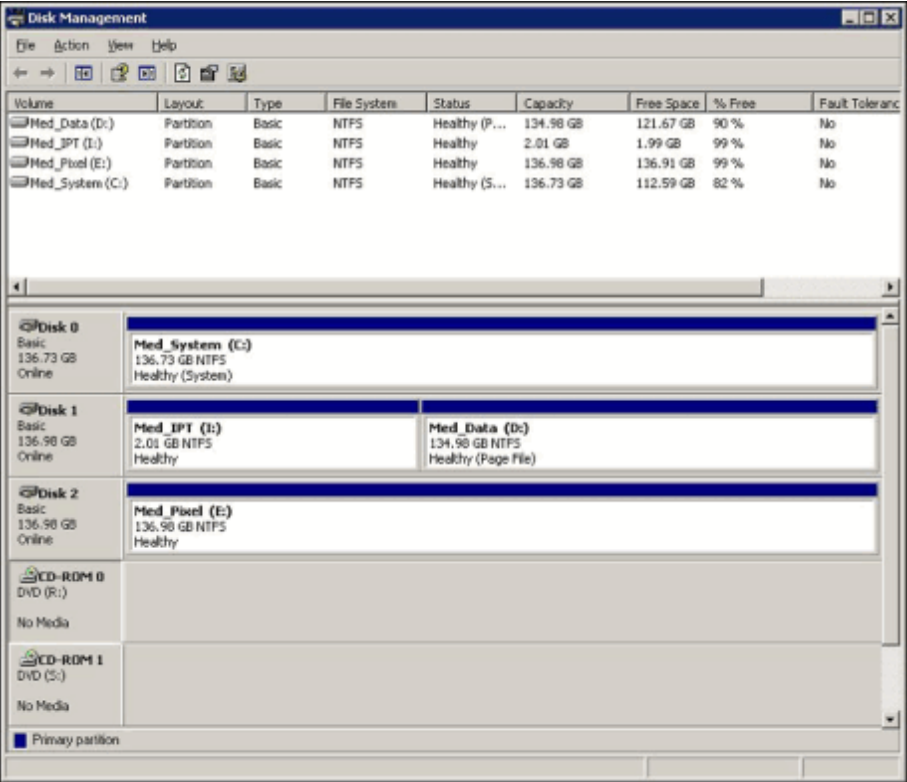


- Perform a reboot.
- Reconfigure Offline devices

5.2.8 Drive letter assignment for hard disks

5.2.8.1 Description:

The system will not start up. The following message pop-up ist displayed.

Fig. 206: HD check message pop-up

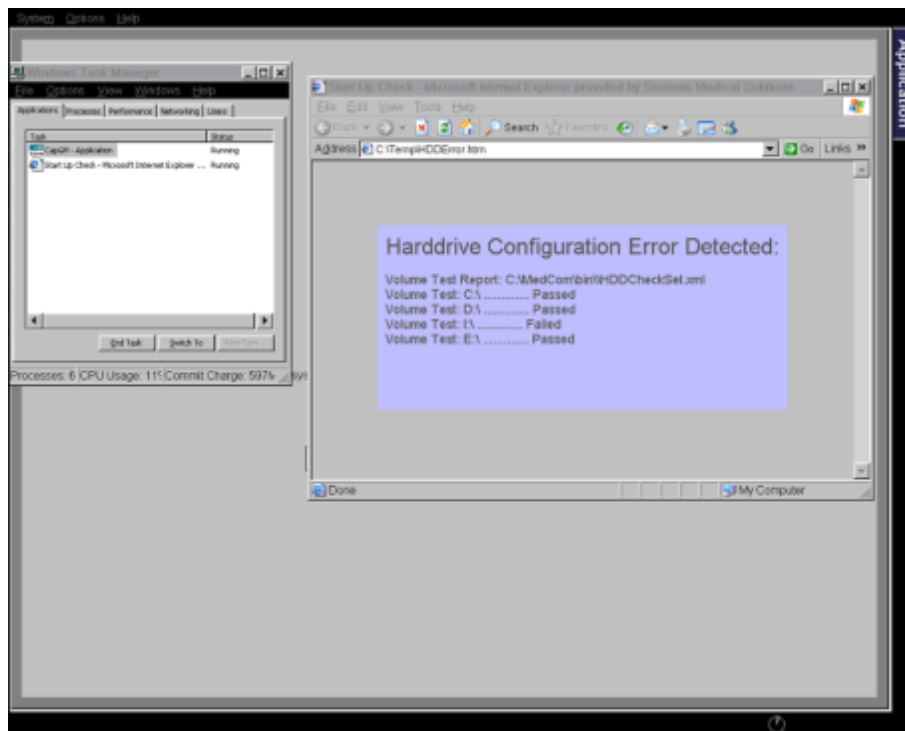


This pop-up shows that the drive letter assignment for the hard disks is not correct.

5.2.8.2 Solution:

The following procedure describes the correction for the drive letter assignment.

- Log in as administrator
- Start the disk manager by selecting **START --> RUN** and entering `diskmgmt.msc`

Fig. 207: RUN pop-up

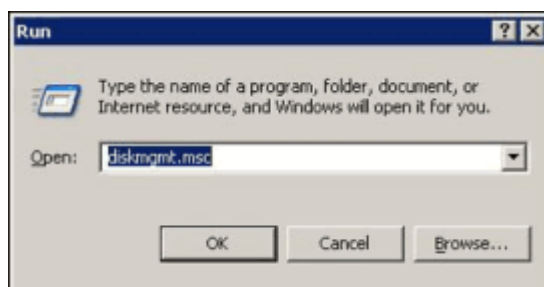

- The disk manager pops up.
- Select the Volume (partition) to be (re-) assigned with a right mouse click, select '**Change drive letter and paths...**' from the pop-up menu, select '**Change**' in the next pop-up and choose the proper drive letter from the drop-down list of the next one. Click '**Yes**' in the 'Confirm' pop-up. (Repeat this sequence for every drive for which you have to change the drive letter).

The following list shows the correct assignment:

- Med_System -> C:
- Med_IPT -> I:

- Med_Data -> D:
- Med_Pixel -> E:

Fig. 208: Drive letters OK - MRC


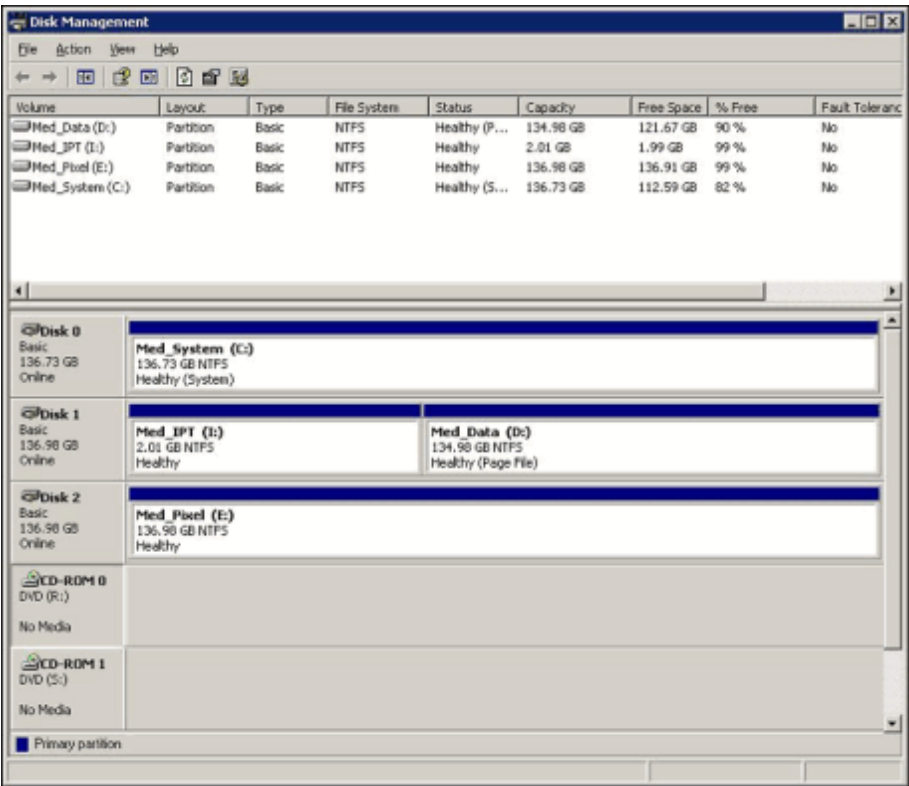


Fig. 209: Drive letters OK - MRSC

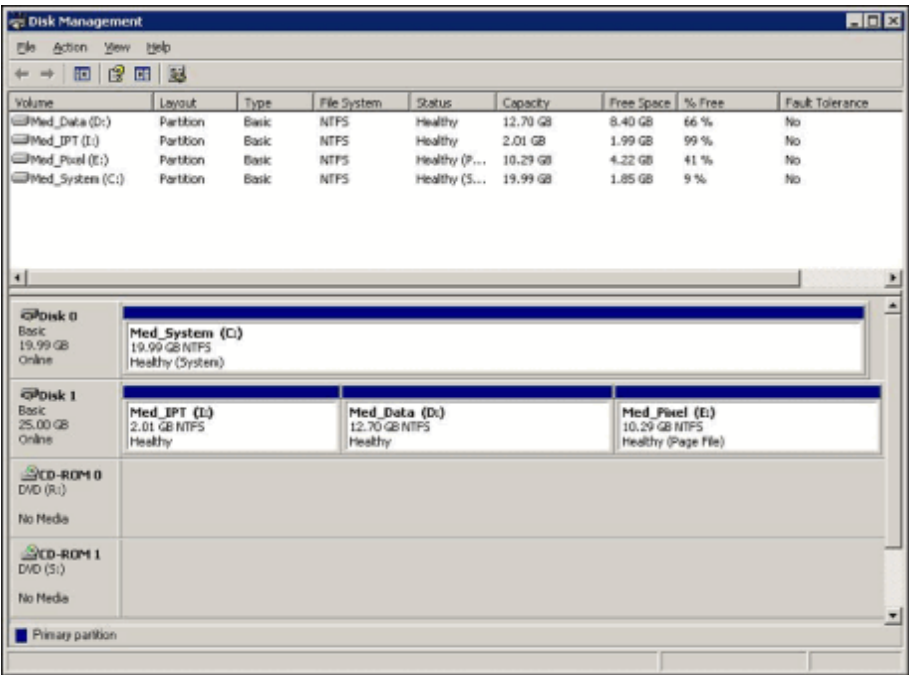

- Perform a reboot.

5.2.9 COBR Backup was interrupted, no System Boot

5.2.9.1 Description:

The backup was interrupted and the system cannot boot from the hard disk. The following message is displayed: "Operating System not found".

5.2.9.2 Solution:

Reset the system to normal boot, do the following:

- Insert the first DVD from the last created backup image into the CD/DVD-ROM drive.
- Wait until the COBR tool window is displayed.
- In the main Backup/Restore pop-up, press Cancel to exit Backup/Restore.
- Confirm the pop-up with Yes to exit.
- The CD/DVD-ROM drive opens, remove the ejected medium from the drive.
- Close the drive and click OK to confirm.
- The system will reboot.
- After the reboot, an autologin as administrator is executed.

After the autologin, additional configuration steps are performed. Wait until they are finished and the system is rebooted.

Do not close or confirm any command windows that are displayed; the system will automatically reboot.

Initial version.

There are no Hazard IDs in this document.

All documents may only be used for rendering services on Siemens Healthcare Products. Any document in electronic form may be printed once. Copy and distribution of electronic documents and hardcopies is prohibited. Offenders will be liable for damages. All other rights are reserved.

healthcare.siemens.com/services

Siemens Healthineers Headquarters
Siemens Healthcare GmbH
Henkestr. 127
91052 Erlangen
Germany
Telephone: +49 9131 84-0
siemens.com/healthineers

Print No.: MR-000.816.48.01.02 | Replaces: n.a.
Doc. Gen. Date: 06.17 | Language: English
© Siemens Healthcare GmbH, 2017

siemens.com/healthineers