

IMAGE DATA SHEET

QA Phantom #1 (cont.)


QA#1/S/10mm/120KV/340mA/1sec . DO NOT use this technique for Vx image series.
Scan at 0 mm

Exam # # of Slice	MTF	MTF 4 slice average	Contrast Scale	Comments	Scan Tech Aprov.
SPECS	---	0.68 to 1.0	110.0 to 130.0	---	---

Box Size : 17 x 17 pixels

Image Acceptance Date: _____ Image Reviewer: _____

QA#1 MTF and CONTRAST SCALE

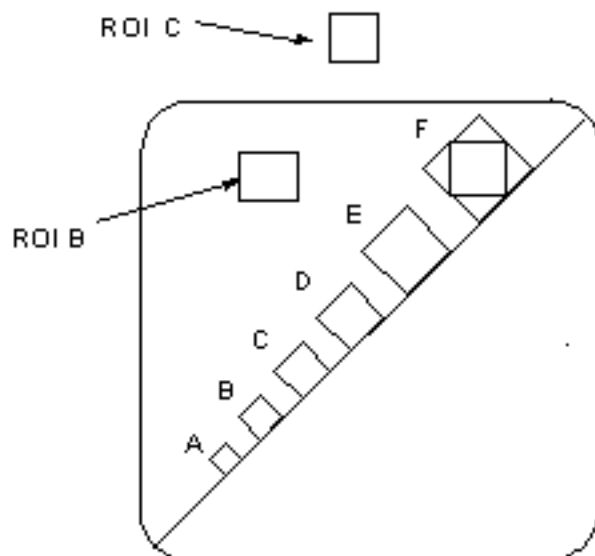

SD_A = Standard deviation of ROI Box over F bar pattern in plexiglass.

SD_{ave} = Average standard deviation of plexiglass and water (i.e., image noise),
 $SD_{ave} = (SD_B + SD_C) / 2$.

Contrast Scale = mean of B - mean of C

Modulation = $(SD_A^2 - SD_{ave}^2)^{1/2}$

MTF = $2.2 \times (\text{Modulation} / \text{Contrast Scale})$