


GE Healthcare

LightSpeed Family

De-Installation Manual

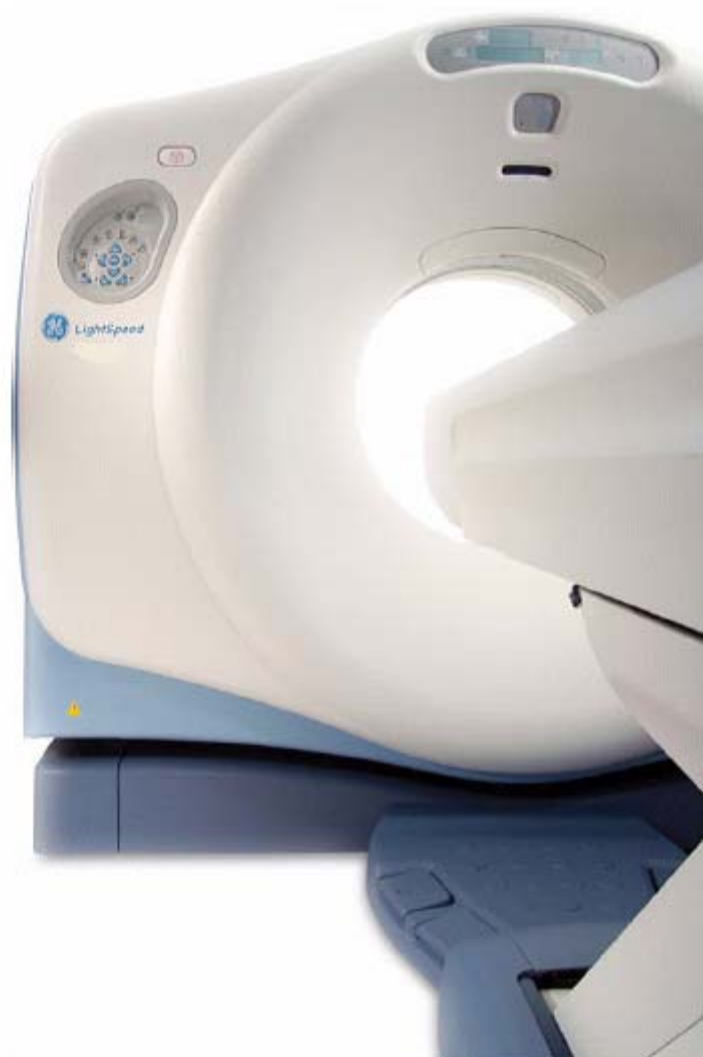

ADVANCED SERVICE DOCUMENTATION. PROPERTY OF GE HEALTHCARE
USE OF THESE MATERIALS LIMITED TO AGENTS AND EMPLOYEES OF
GE HEALTHCARE OR OTHER PARTIES EXPRESSLY LICENSED BY
GE HEALTHCARE. UNLICENSED USE IS STRICTLY PROHIBITED.



5272259-1EN

Rev 1

© 2008, General Electric Company,
All Rights Reserved

IMPORTANT PRECAUTIONS**LANGUAGE****ПРЕДУПРЕЖДЕНИЕ**

(BG)

- ТОВА УПЪТВАНЕ ЗА РАБОТА Е НАЛИЧНО САМО НА АНГЛИЙСКИ ЕЗИК.
- АКО ДОСТАВЧИКЪТ НА УСЛУГАТА НА КЛИЕНТА ИЗИСКА ЕЗИК, РАЗЛИЧЕН ОТ АНГЛИЙСКИ, ЗАДЪЛЖЕНИЕ НА КЛИЕНТА Е ДА ОСИГУРИ ПРЕВОД.
- НЕ ИЗПОЛЗВАЙТЕ ОБОРУДВАНЕТО ПРЕДИ ДА СТЕ СЕ КОНСУЛТИРАЛИ И РАЗБРАЛИ УПЪТВАНЕТО ЗА РАБОТА.
- НЕСПАЗВАНЕТО НА ТОВА ПРЕДУПРЕЖДЕНИЕ МОЖЕ ДА ДОВЕДЕ ДО НАРАНЯВАНЕ НА ДОСТАВЧИКА НА УСЛУГАТА, ОПЕРАТОРА ИЛИ ПАЦИЕНТ В РЕЗУЛТАТ НА ТОКОВ УДАР ИЛИ МЕХАНИЧНА ИЛИ ДРУГА ОПАСНОСТ.

警告

(ZH-CN)

- 本维修手册仅提供英文版本。
- 如果维修服务提供商需要非英文版本，客户需自行提供翻译服务。
- 未详细阅读和完全理解本维修手册之前，不得进行维修。
- 忽略本警告可能对维修人员，操作员或患者造成触电、机械伤害或其他形式的伤害。

VÝSTRAHA

(CS)

- TENTO PROVOZNÍ NÁVOD EXISTUJE POUZE V ANGLICKÉM JAZYCE.
- V PŘÍPADĚ, ŽE EXTERNÍ SLUŽBA ZÁKAZNÍKŮM POTŘEBUJE NÁVOD V JINÉM JAZYCE, JE ZAJIŠTĚNÍ PŘEKladU DO ODPOVÍDAJÍCÍHO JAZYKA ÚKOLEM ZÁKAZNÍKA.
- NESNAŽTE SE O ÚDRŽBU TOHOTO ZAŘÍZENÍ, ANIŽ BYSTE SI PŘEČETLI TENTO PROVOZNÍ NÁVOD A POCHOPILI JEHO OBSAH.
- V PŘÍPADĚ NEDODRŽOVÁNÍ TÉTO VÝSTRAHY MŮŽE DOJÍT K PORANĚNÍ PRACOVNÍKA PRODEJNÍHO SERVISU, OBSLUŽNÉHO PERSONÁLU NEBO PACIENTŮ VLIVEM ELEKTRICKÉHO PROUDU, RESPEKTIVE VLIVEM MECHANICKÝCH ČI JINÝCH RIZIK.

ADVARSEL

(DA)

- DENNE SERVICEMANUAL FINDES KUN PÅ ENGELSK.
- HVIS EN KUNDES TEKNIKER HAR BRUG FOR ET ANDET SPROG END ENGELSK, ER DET KUNDENS ANSVAR AT SØRGE FOR OVERSÆTTELSE.
- FORSØG IKKE AT SERVICERE UDSTYRET MEDMINDRE DENNE SERVICEMANUAL HAR VÆRET KONSULTERET OG ER FORSTÅET.
- MANGLENDE OVERHOLDELSE AF DENNE ADVARSEL KAN MEDFØRE SKADE PÅ GRUND AF ELEKTRISK, MEKANISK ELLER ANDEN FARE FOR TEKNIKEREN, OPERATØREN ELLER PATIENTEN.

WAARSCHUWING

(NL)

- Deze onderhoudshandleiding is enkel in het Engels verkrijgbaar.
- Als het onderhoudspersoneel een andere taal vereist, dan is de klant verantwoordelijk voor de vertaling ervan.
- Probeer de apparatuur niet te onderhouden voordat deze onderhoudshandleiding werd geraadpleegd en begrepen is.
- Indien deze waarschuwing niet wordt opgevolgd, zou het onderhoudspersoneel, de operator of een patiënt gewond kunnen raken als gevolg van een elektrische schok, mechanische of andere gevaren.

WARNING

(EN)

- This Service Manual is available in English only.
- If a customer's service provider requires a language other than English, it is the customer's responsibility to provide translation services.
- Do not attempt to service the equipment unless this service manual has been consulted and is understood.
- Failure to heed this warning may result in injury to the service provider, operator, or patient from electric shock, or from mechanical or other hazards.

HOIATUS

(ET)

- Käesolev teenindusjuhend on saadaval ainult inglise keeles.
- Kui klienditeeninduse osutaja nõuab juhendit inglise keelest erinevas keeles, vastutab klient tõlketeenuse osutamise eest.
- Ärge üritage seadmeid teenindada enne eelnevalt käesoleva teenindusjuhendiga tutvumist ja sellest aru saamist.
- Käesoleva hoiatuse eiramine võib põhjustada teenuseosutaja, operaatori või patsiendi vigastamist elektrilöögi, mehaanilise või muu ohu tagajärjel.

VAROITUS

(FI)

- Tämä huolto-ohje on saatavilla vain englanniksi.
- Jos asiakkaan huoltohenkilöstö vaatii muuta kuin englanninkielistä materiaalia, tarvittavan käännöksen hankkiminen on asiakkaan vastuulla.
- Älä yritä korjata laitteistoa ennen kuin olet varmasti lukenut ja ymmärtänyt tämän huolto-ohjeen.
- Mikäli tätä varoitusta ei noudateta, seurauksena voi olla huoltohenkilöstön, laitteiston käyttäjän tai potilaan vahingoittuminen sähköiskun, mekaanisen vian tai muun vaaratilanteen vuoksi.

ATTENTION

(FR)

- Ce manuel de service n'est disponible qu'en anglais.
- Si le technicien du client a besoin de ce manuel dans une autre langue que l'anglais, c'est au client qu'il incombe de le faire traduire.
- Ne pas tenter d'intervenir sur les équipements tant que le manuel service n'a pas été consulté et compris.
- Le non-respect de cet avertissement peut entraîner chez le technicien, l'opérateur ou le patient des blessures dues à des dangers électriques, mécaniques ou autres.

WARNUNG

(DE)

- Diese Serviceanleitung existiert nur in Englischer Sprache.
- Falls ein fremder Kundendienst eine andere Sprache benötigt, ist es Aufgabe des Kunden für eine Entsprechende Übersetzung zu sorgen.
- Versuchen Sie nicht diese Anlage zu warten, ohne diese Serviceanleitung gelesen und verstanden zu haben.
- Wird diese Warnung nicht beachtet, so kann es zu Verletzungen des Kundendiensttechnikers, des Bedieners oder des Patienten durch Stromschläge, Mechanische oder Sonstige Gefahren kommen.

ΠΡΟΕΙΔΟΠΟΙΗΣΗ
(EL)

- Το παρόν εγχειρίδιο σέρβις διατίθεται στα αγγλικά μόνο.
- Εάν το άτομο παροχής σέρβις ενός πελάτη απαιτεί το παρόν εγχειρίδιο σε γλώσσα εκτός των αγγλικών, αποτελεί ευθύνη του πελάτη να παρέχει υπηρεσίες μετάφρασης.
- Μην επιχειρήσετε την εκτέλεση εργασιών σέρβις στον εξοπλισμό εκτός εάν έχετε συμβουλευτεί και έχετε κατανοήσει το παρόν εγχειρίδιο σέρβις.
- Εάν δε λάβετε υπόψη την προειδοποίηση αυτή, ενδέχεται να προκληθεί τραυματισμός στο άτομο παροχής σέρβις, στο χειριστή ή στον ασθενή από ηλεκτροπληξία, μηχανικούς ή άλλους κινδύνους.

FIGYELMEZTETÉS
(HU)

- Ezen karbantartási kézikönyv kizárólag angol nyelven érhető el.
- Ha a vevő szolgáltatója angoltól eltérő nyelvre tart igényt, akkor a vevő felelőssége a fordítás elkészíttetése.
- Ne próbálja elkezdni használni a berendezést, amíg a karbantartási kézikönyvben leírtakat nem értelmezték.
- Ezen figyelmeztetés figyelmen kívül hagyása a szolgáltató, működtető vagy a beteg áramütés, mechanikai vagy egyéb veszélyhelyzet miatti sérülését eredményezheti.

AÐVÖRUN
(IS)

- Þessi þjónustuhandbók er eingöngu fáanleg á ensku.
- Ef að þjónustuveitandi viðskiptamanns þarfnast annas tungumáls en ensku, er það skylda viðskiptamanns að skaffa tungumálþjónustu.
- Reynið ekki að afgreiða tækið nema að þessi þjónustuhandbók hefur verið skoðuð og skilin.
- Brot á þessari aðvörun getur leitt til meiðsla á þjónustuveitanda, stjórnanda eða sjúklings frá raflosti, vélrænu eða öðrum áhættum.

AVVERTENZA
(IT)

- Il presente manuale di manutenzione è disponibile soltanto in inglese.
- Se un addetto alla manutenzione richiede il manuale in una lingua diversa, il cliente è tenuto a provvedere direttamente alla traduzione.
- Si proceda alla manutenzione dell'apparecchiatura solo dopo aver consultato il presente manuale ed averne compreso il contenuto.
- Il non rispetto della presente avvertenza potrebbe far compiere operazioni da cui derivino lesioni all'addetto, alla manutenzione, all'utilizzatore ed al paziente per folgorazione elettrica, per urti meccanici od altri rischi.

警告
(JA)

- このサービスマニュアルには英語版しかありません。
- サービスを担当される業者が英語以外の言語を要求される場合、翻訳作業はその業者の責任で行うものとさせていただきます。
- このサービスマニュアルを熟読し理解せずに、装置のサービスを行わないでください。
- この警告に従わない場合、サービスを担当される方、操作員あるいは患者さんが、感電や機械的又はその他の危険により負傷する可能性があります。

경고
(KO)

- 본 서비스 지침서는 영어로만 이용하실 수 있습니다.
- 고객의 서비스 제공자가 영어 이외의 언어를 요구할 경우, 번역 서비스를 제공하는 것은 고객의 책임입니다.
- 본 서비스 지침서를 참고했고 이해하지 않는 한은 해당 장비를 수리하려고 시도하지 마십시오.
- 이 경고에 유의하지 않으면 전기 쇼크, 기계상의 혹은 다른 위험으로부터 서비스 제공자, 운전자 혹은 환자에게 위해를 가할 수 있습니다.

BRĪDINĀJUMS

(LV)

- Šī apkalpes rokasgrāmata ir pieejama tikai angļu valodā.
- Ja klienta apkalpes sniedzējam nepieciešama informācija citā valodā, nevis angļu, klienta pienākums ir nodrošināt tulkošanu.
- Neveiciet aprikojuma apkalpi bez apkalpes rokasgrāmatas izlasīšanas un saprašanas.
- Šī brīdinājuma neievērošana var radīt elektriskās strāvas trieciena, mehānisku vai citu risku izraisītu traumu apkalpes sniedzējam, operatoram vai pacientam.

ISPĒJIMAS

(LT)

- Šīs ekspluatavimo vadovas yra prieinamas tik anglų kalba.
- Jei kliento paslaugų tiekėjas reikalauja vadovo kita kalba – ne anglų, numatyti vertimo paslaugas yra kliento atsakomybė.
- Nemėginkite atlikti įrangos techninės priežiūros, nebent atsižvelgėte į šį ekspluatavimo vadovą ir jį supratote.
- Jei neatkreipsite dėmesio į šį perspėjimą, galimi sužalojimai dėl elektros šoko, mechaninių ar kitų pavojų paslaugų tiekėjui, operatoriui ar pacientui.

ADVARSEL

(NO)

- Denne servicehåndboken finnes bare på engelsk.
- Hvis kundens serviceleverandør trenger et annet språk, er det kundens ansvar å sørge for oversettelse.
- Ikke forsøk å reparere utstyret uten at denne servicehåndboken er lest og forstått.
- Manglende hensyn til denne advarselen kan føre til at serviceleverandøren, operatøren eller pasienten skades på grunn av elektrisk støt, mekaniske eller andre farer.

OSTRZEŻENIE

(PL)

- Niniejszy podręcznik serwisowy dostępny jest jedynie w języku angielskim.
- Jeśli dostawca usług klienta wymaga języka innego niż angielski, zapewnienie usługi tłumaczenia jest obowiązkiem klienta.
- Nie próbować serwisować wyposażenia bez zapoznania się i zrozumienia niniejszego podręcznika serwisowego.
- Niezastosowanie się do tego ostrzeżenia może spowodować urazy dostawcy usług, operatora lub pacjenta w wyniku porażenia elektrycznego, zagrożenia mechanicznego bądź innego.

ATENÇÃO

(PT)

- Este manual de assistência técnica só se encontra disponível em inglês.
- Se qualquer outro serviço de assistência técnica solicitar estes manuais noutro idioma, é da responsabilidade do cliente fornecer os serviços de tradução.
- Não tente consertar o equipamento sem ter consultado e compreendido este manual de assistência técnica
- O não cumprimento deste aviso pode pôr em perigo a segurança do técnico, do operador ou do paciente devido a choques elétricos, mecânicos ou outros.

ATENȚIE

(RO)

- Acest manual de service este disponibil numai în limba engleză.
- Dacă un furnizor de servicii pentru clienți necesită o altă limbă decât cea engleză, este de datoria clientului să furnizeze o traducere.
- Nu încercați să reparați echipamentul decât ulterior consultării și înțelegerii acestui manual de service.
- Ignorarea acestui avertisment ar putea duce la rănirea depanatorului, operatorului sau pacientului în urma pericolelor de electrocutare, mecanice sau de altă natură.

ОСТОРОЖНО!
(RU)

- Данное руководство по обслуживанию предлагается только на английском языке.
- Если сервисному персоналу клиента необходимо руководство не на английском, а на каком-то другом языке, клиенту следует самостоятельно обеспечить перевод.
- Перед обслуживанием оборудования обязательно обратитесь к данному руководству и поймите изложенные в нем сведения.
- Несоблюдение требований данного предупреждения может привести к тому, что специалист по обслуживанию, оператор или пациент получат удар электрическим током, механическую травму или другое повреждение.

UPOZORNENIE
(SK)

- Tento návod na obsluhu je k dispozícii len v angličtine.
- Ak zákazníkovi poskytovateľ služieb vyžaduje iný jazyk ako angličtinu, poskytnutie prekladateľských služieb je zodpovednosťou zákazníka.
- Nepokúšajte sa o obsluhu zariadenia skôr, ako si neprečítate návod na obsluhu a neporozumiete mu.
- Zanedbanie tohto upozornenia môže vyústiť do zranenia poskytovateľa služieb, obsluhujúcej osoby alebo pacienta elektrickým prúdom, do mechanického alebo iného nebezpečenstva.

ATENCION
(ES)

- Este manual de servicio sólo existe en inglés.
- Si el encargado de mantenimiento de un cliente necesita un idioma que no sea el inglés, el cliente deberá encargarse de la traducción del manual.
- No se deberá dar servicio técnico al equipo, sin haber consultado y comprendido este manual de servicio.
- La no observancia del presente aviso puede dar lugar a que el proveedor de servicios, el operador o el paciente sufran lesiones provocadas por causas eléctricas, mecánicas o de otra naturaleza.

VARNING
(SV)

- Den här servicehandboken finns bara tillgänglig på engelska.
- Om en kunds servicetekniker har behov av ett annat språk än engelska ansvarar kunden för att tillhandahålla översättningstjänster.
- Försök inte utföra service på utrustningen om du inte har läst och förstår den här servicehandboken.
- Om du inte tar hänsyn till den här varningen kan det resultera i skador på serviceteknikern, operatören eller patienten till följd av elektriska stötar, mekaniska faror eller andra faror.

DIKKAT
(TR)

- Bu servis kılavuzunun sadece ingilizcesi mevcuttur.
- Eğer müşteri teknisyeni bu kılavuzu ingilizce dışında bir başka lisandan talep ederse, bunu tercüme ettirmek müşteriye düşer.
- Servis kılavuzunu okuyup anlamadan ekipmanlara müdahale etmeyiniz.
- Bu uyarıya uyulmaması, elektrik, mekanik veya diğer tehlikelerden dolayı teknisyen, operatör veya hastanın yaralanmasına yol açabilir.

CERTIFIED ELECTRICAL CONTRACTOR STATEMENT

All licensed electrical contractors shall perform all electrical de-installation work required to remove electrical feeds from the Power Distribution Unit and the A1 disconnect. The products involved are highly sophisticated, and special engineering competence is required. All de-installation electrical work must comply with the requirements of the applicable electrical codes. GE Healthcare (GEH) shall only utilize trained personnel (i.e., General Electric field engineers, personnel of third-party service companies with equivalent training, or licensed electricians) to perform the de-installation work.

IMPORTANT...X-RAY PROTECTION

X-ray equipment if not properly used may cause injury. Although this apparatus incorporates a high degree of protection against x-radiation other than the useful beam, no practical design of equipment can provide complete protection. Nor can any practical design compel the operator to take adequate precautions to prevent the possibility of any persons carelessly exposing themselves or others to radiation.

It is important that anyone having anything to do with x-radiation be properly trained and fully acquainted with the fundamentals of radiation, completed a course in Radiation Protection, and take adequate steps to protect themselves against injury. Various protective materials and devices are available. It is urged that such materials or devices be used.

The equipment is sold with the understanding that the General Electric Company, GE Healthcare, its agents, and representatives have no responsibility for injury or damage which may result from improper use of the equipment.

OMISSIONS & ERRORS

Customers, please contact your GE Sales or Service representatives.

GE personnel, please use the GE Healthcare iTrak Process to report all omissions, errors, and defects in this publication.

REVISION HISTORY

Rev	Date	Author	Reason For Change
1	May 6, 2008	J. Neumiller	Initial version.

TABLE OF CONTENTS

REVISION HISTORY.....	8
TABLE OF CONTENTS.....	9
CHAPTER 1- INTRODUCTION	11
1-1 Skill Set Required.....	11
1-2 De-installation Leader Requirements.....	11
1-3 De-Installation Safety.....	12
1-4 General De-Installation Best Practices.....	12
1-5 Required Common Tools and Supplies.....	12
Wrenches	12
Sockets and Extensions	12
Screw Drivers	13
Power Tools.....	13
Electrical Tools	13
Hand Tools	13
System Cleaners	13
Personal Protective Equipment.....	13
CHAPTER 2 – BEFORE YOU START DE-INSTALLING EQUIPMENT.....	14
2-1 Additional Documents Required	14
2-1-1 Returning Used Systems	14
2-2 Equipment Inventory.....	14
2-3 Floor Protection	14
2-4 Equipment Handling	14
2-4-1 Lifting Requirements	14
2-4-2 Door Openings.....	15
2-4-3 Elevator Requirements	15
2-5 HIPAA	15
2-5-1 Patient Privacy Notice.....	15
2-6 Customer Task List for De-installation.....	15
2-7 De-Installation Prerequisites.....	16
2-7-1 Tube Cooling.....	16
2-7-2 Decontamination	16
2-8 Console Clean-up and Power-down.....	16
CHAPTER 3- PREPARE CONSOLE FOR SHIPMENT	17
3-1 Field Engineer De-installation Tasks	17
3-2 System Save State *.....	17
3-3 Shut Down Applications.....	17
3-4 Remove Customer (Hospital/Clinic) Information / Reconfigure Operator Console	18
3-5 Remove Customer Items	19
3-5-1 Remove Customer's MODs and DVDs.....	19
3-6 Power-down Console.....	19
CHAPTER 4- DE-INSTALLATION POWER DOWN & LOTO.....	20
4-1 A-1 Disconnect Lock Out Procedure	20
CHAPTER 5- MECHANICAL DE-INSTALLATION OVERVIEW.....	21
CHAPTER 6- DE-INSTALL THE OPERATOR CONSOLE.....	22
6-1 Remove & Pack Monitors, SCIM Assembly and SCSI Tower.....	22
6-2 Disassemble & Pack Keyboard Table Top Accessories	22
6-3 VCT ONLY- Remove Hard Drives	22
6-4 Remove Console Rear Cable Connections.....	22
6-5 Remove all Disconnected System Cables.....	22

6-6 Pack console for shipping.....	22
7-1 De-Install Power Distribution Unit.....	24
7-2 De-Install Uninterrupted Power Supply (UPS).....	24
CHAPTER 8- DE-INSTALLATION OF THE TABLE	26
CHAPTER 9- DE-INSTALLATION OF THE GANTRY	27
CHAPTER 10- SYSTEM OPTIONS	29
10-1 Injectors	29
10-1-1 Floor Mounted Injectors	29
10-2 Monitors	29
10-3 RPM D4D - Respiratory Gating	29
CHAPTER 11- CABLES, ACCESSORIES AND FINAL PROCEDURES	31
CHAPTER 11- CABLES, ACCESSORIES AND FINAL PROCEDURES	31
APPENDIX A	32
Customer Items.....	33

Chapter 1- Introduction

1-1 Skill Set Required

All GE Healthcare (GEH) representatives must meet the minimum requirements to perform de-installation work. GEH shall only utilize trained personnel (i.e., GEH field engineers, personnel of third-party service companies with equivalent training, and licensed electricians) to perform system de-installations.

Note: Before you start the de-installation, make sure the site preparation complies with conditions and instructions found in this manual.

1-2 De-installation Leader Requirements

To lead a de-installation one must be trained as a Field Engineer (FE) or Supplier to install LightSpeed, BrightSpeed, RT, or equivalent systems, and have experience with one of the following:

- A full system installation as a leader.
- A de-installation as a helper.

Required training courses include:

- Environmental Health and Safety (EHS)
- Hazardous Communication
- Eye Protection, and Personal Protection Equipment (PPE)
- Blood Borne Pathogens
- Lock-out Tag-out procedure
- Professional Communication

A FE or Supplier may meet the requirements through:

- My Learning training verification and/or on-demand from supplier.
- Each Supplier is responsible for providing the required training.
- Each Supplier's team member shall have all required PPE as listed, while on site.

1-3 De-Installation Safety

Before performing a de-installation, review system safety procedures in the service manual. During the de-installation one could come in contact with several types of hazards that can cause bodily injuries. Follow all safety requirements when working around and in contaminated areas.

WARNING

BIOHAZARDS COULD BE PRESENT WHEN REMOVING COMPONENTS AND COVERS, WHICH HAVE BEEN EXPOSED TO PATIENT CONTACT. FOLLOW YOUR ORGANIZATION'S BLOOD BOURNE PATHOGENE PROCEDURES, AND POLICIES. FOLLOW ALL PERSONAL PROTECTIVE EQUIPMENT REQUIREMENTS.

Bio-hazardous areas can be found around the cradle surface, the table covers, the gantry scan window and base covers. Hazardous items such as broken glass, needles, and bodily fluids are often found behind and under the table foot switch and gantry base covers.

Biohazards involve infectious biological agents present in human blood, bodily fluids and blood products such as:

- Hepatitis B Virus (HBV)
- Human Immune-deficiency Virus (HIV)
- Tuberculosis bacterium

1-4 General De-Installation Best Practices

- Do not accept equipment for removal that **has not been cleaned** by the user. If the area is contaminated, have the staff clean the area before working in the area.
- Always wear gloves during the de-installation process.
- Exercise caution when removing components that have contact with body fluids such as whole blood.
- Wash your hands before going from a contaminated area to a non-contaminated area, or when gloves are removed or changed.
- Keep your hands away from your face.
- Cover all superficial cuts and wounds before starting work.
- Disinfect tools and other items that have been near any part of the instrument sample path or waste area with 10% v/v bleach after washing with warm soapy water.

1-5 Required Common Tools and Supplies

The following tools and supplies are required for a complete de-installation.

Wrenches

Standard and metric combination wrench sets
Standard and metric hex key (Allen wrench) sets

Sockets and Extensions

3/8" drive metric and standard socket set
1", 1-1/8", 1-1/4", 1-1/2", 1-5/8" sockets
1/2" drive ratchet wrenches
3/4" deep well socket
Metric hex bit set 1/4" or 3/8" drive, including 10mm and 14mm hex bit

3/8" drive universal joint
21mm socket (optional)

Screw Drivers

Phillips screwdriver set (small, medium, and large)
Straight blade screwdriver set (small, medium, and large)
Pozi-drive 0 screwdriver
Pozi-drive #1 screwdriver

Power Tools

3/8" or 1/2" drill, cordless or electric
Reciprocating saw (Sawzall or equivalent) and assorted blades or
5" grinder with metal cutting wheel
25' extension power cords

Electrical Tools

DVM
Monitor Static Kit

Hand Tools

Ball-peen hammer (1lb or 2lb)
Diagonal cutting pliers, small and large
Tongue & groove pliers (large)
Large pry bar
Diagonal cutting pliers, large (to cut 1/0 ground)

System Cleaners

Purchase Locally:
409 Antibacterial Cleaner or equivalent
Hand towels
Hand sanitizer

Personal Protective Equipment

Safety shoes*
LOTO kit* -- MUST have **approved GE** LOTO kit
Safety glasses*
Leather gloves*
Hearing protection*
Knee pads
6' or 8' step ladders
Rope or strap with a load rating of 300#

* Required items

Chapter 2 – Before You Start De-Installing Equipment

2-1 Additional Documents Required

- Installation manual
- Pre-installation manual

Locate the pre-installation manual and installation manual **shipped with the system**. If the paper copy is not available, both manuals can be located on the white service CD ROM shipped with the system. Check the service cabinet for these items.

Review the delivery and storage sections in the pre-installation manual. Review the installation manual assembly instructions **prior** to disassembly of the system.

2-1-1 Returning Used Systems

Kits are available from [UMI Website](#). Listed under GoldSeal returned packaging

CT Kit

- # 10 for H.S.A. Systems
- # 14 for LightSpeed and BrightSpeed Families

If you are not returning the system to GE Healthcare for recycling, be sure to comply with local and national regulations for proper disposal and recycling.

2-2 Equipment Inventory

The CT system ships in multiple boxes, on pallets and dollies, with approximately 26 units. Begin by locating the De-installation kit and system dollies. An example inventory sheet is located in Appendix A.

2-3 Floor Protection

The de-installation team is required to provide floor protection in the room where the equipment is being de-installed. The movers are responsible for floor protection from the scan room to the loading area of their truck.

2-4 Equipment Handling

2-4-1 Lifting Requirements

The Project Manager of Installation (PMI) should contact engineering for all special lifting requirements.

The following items must **never** be dropped:

- Gantry
- Console
- Table
- PDU

A drop from a height greater than one half inch ($\frac{1}{2}$ " or 12mm) will induce structural damage to the frame or other major components. This damage will not always be obvious or clearly visible.

All system components must remain upright at all times, and must not be tipped.

Unauthorized gantry lifting will cause gantry-bearing damage:

- The gantry cannot be hoisted using the shipping dollies.
- The gantry cannot be lifted with a forklift.
- All gantry lifting requires engineering approval for each occurrence.

2-4-2 Door Openings

Door openings for moving equipment through building must be 42" X 82" (1067mm X 2083mm), and a 8' (2439mm) corridor width is recommended.

2-4-3 Elevator Requirements

Refer to the Pre-Installation Manual to determine equipment specifications and weights.

2-5 HIPAA

In the USA, the *Health and Human Services* (HHS) department has established a set of standards entitled *Health Insurance Portability and Accountability Act* (HIPAA). A key component of these standards are the laws and ensuing regulations, which established guidelines for individual patient management and the access to healthcare information. Included is the provision of a legal firewall to prevent inappropriate, unauthorized disclosure of patient information. This applies to access of confidential electronic health information found on the console computer and hard drives, PACC systems, and Advantage Workstations.

2-5-1 Patient Privacy Notice

Laws and regulations are established and globally enforced in all countries for patient privacy management and unauthorized access to patient healthcare information.

In an effort to prevent inadvertent disclosures of patient information during a system de-installation, GE Healthcare and/or their agents shall remove all patient data, patient images, and image raw data prior to system final power down at de-installation. Procedures in this manual represent the minimum tasks required to meet global patient privacy standards, acts and regulations.

If your country, state, region or local laws regarding patient healthcare and privacy standards are more stringent than stated in this manual, always follow the requirements for your area.

Included Sub-systems

- Console computers and hard drives
- PACC systems and storage systems
- Advantage Workstations and storage systems

2-6 Customer Task List for De-installation

The PMI will notify the customer that prior to the de-installation of their system, the following tasks should be completed:

1. Remove all images from the system.
2. Save images on MOD or media of choice.
3. Empty all personal contents from console drawer.
4. Store all hospital records and data in a secure area.

2-7 De-Installation Prerequisites

Confirm the conditions listed here have been reviewed and actions taken.

2-7-1 Tube Cooling

The tube must be cool before de-installation of the gantry. Wait one hour after the last system scan before performing work inside the gantry. If the tube is hot continue to wait until the tube has cooled.

2-7-2 Decontamination

Verify the area around the system work area is safe, and free of biohazards. If not, follow these procedures:

Have the staff clean the area before working.

Always wear gloves during the de-installation process.

Disinfect tools and other items that have been used in the de-installation process area with 100% v/v bleach, after washing with warm soapy water.

Leather gloves should be worn when working around the table footswitch areas to prevent skin puncture due to sharp edges and debris.

Discard leather gloves after in contaminated areas.

2-8 Console Clean-up and Power-down

The PMI and/or customer's approval must be granted to start the de-installation process. If it is necessary to stop de-installation, once started, the time to return the system to normal operation can be timely.

Confirm the system you plan to de-install is correct as shown on the "notice to proceed" paperwork. If you or the customer has questions, escalate this to the PMI.

The customer should have all images removed from the system, and saved imaged on MOD from the SCSI tower. If customer patient information is still present on the console, ask the customer to remove all patient information.

Chapter 3- Prepare Console for Shipment

This outlines the tasks for the field engineer.

3-1 Field Engineer De-installation Tasks

The field engineer will perform these tasks at a minimum prior to system de-installation:

1. Remove customer image data
2. Remove customer hospital information
3. Regenerate image database
4. Confirm patient and customer information in **not** present

Note: Do not return MOD or DVD disk with patient images or data to GE Healthcare or to the recycle center.

3-2 System Save State *

*Perform this procedure only when the system will be transferred and re-installed in a new location.

The system state backup saves vital calibration and scanner files. Follow the save and restore procedures throughout this document precisely or valuable information may be lost.

Note: This process takes approximately 10 minutes, but in some cases it may take up to 30 minutes.

1. Insert the System Save State disc into the drive.
2. Select: *Service* icon to access the CSD (Common Service Desktop).
3. Select: [Utilities].
4. Select: [System State]. The System State screen appears.
5. Select [All] to select all the calibrations, characterizations, etc.
6. Select [Save]. The Save System State box appears.
7. Select [Yes].
8. Verify that the "Save" of System State was successful. If not, correct any errors and re-save the system state.

A message at the end of the Save should state: *Save/Restore System State: Completed Successfully*

9. Pack disk with console materials.


3-3 Shut Down Applications

1. If Applications is currently running, you must shutdown system applications.
2. Click on the SERVICE DESKTOP icon.
3. On the desktop toolbar, select UTILITIES icon.
4. Click APPLICATIONS SHUTDOWN (to bring down applications only).

3-4 Remove Customer (Hospital/Clinic) Information / Reconfigure Operator Console

1. Shut down applications from the Service Desktop.
2. In an **xterm** window, log in as root:
 - a. Type: `[su -]` and press [Enter]
 - b. Type the `root password [xxxx]` and press [Enter]
3. Launch the **Install** utility
 - a. Type: `[reconfig]` and select [ENTER] at the prompt.
 - b. The Operator Console displays the Install Utility Window
4. Select [Config] and the Operator Console displays the CT Software Configuration Screen-System Settings Screen. See **Illustration 3-1**.
 - a. Remove the customer name, replace with [GEMS].
 - b. Remove the Service ID, replace with [GEMS].
 - c. Select regenerate data base, Enter [Yes].

Illustration 3-1: CT Software Configuration Screen

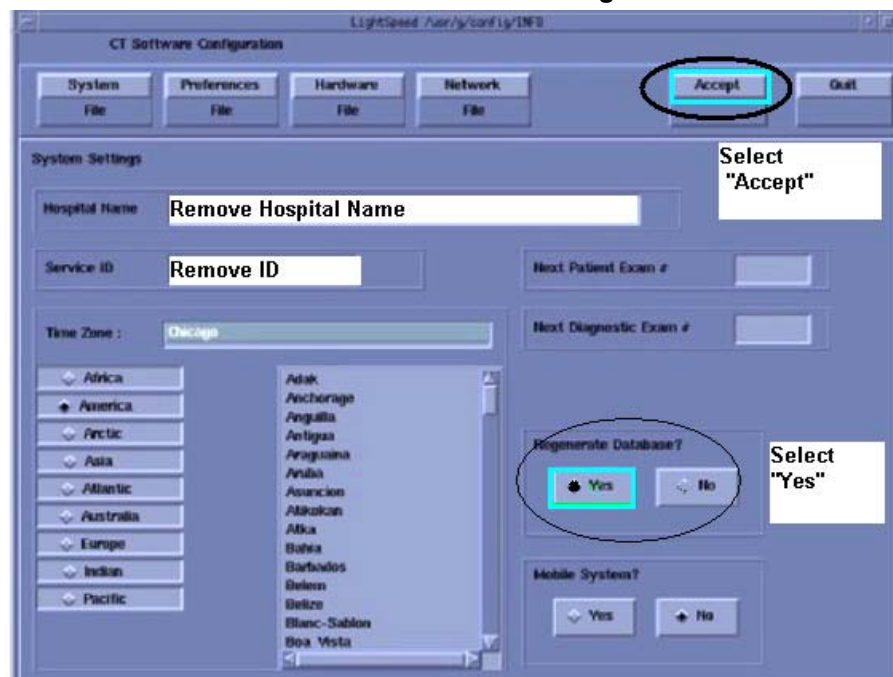

5. Select [ACCEPT]. Be patient this takes a few minutes.
6. Select [YES] to reboot. The system will automatically log in as "ctuser" after the reboot.
7. Select [OK] on the **Autostart Disabled** popup message.

8. Open a Shell window.
 - a. Type: [ST] select [Enter].
9. When the system is at Applications, select [Image Works].
 - a. Confirm that customer image data is not present.

3-5 Remove Customer Items

Locate an empty box and label it "CUSTOMER MATERIALS". Place all remaining customer materials in this box such as image disk, and hospital accessories.

WARNING: Do not return customer materials with the system.

3-5-1 Remove Customer's MODs and DVDs

Remove all MODs and DVDs from the SCSI tower drives.

- Press the eject button on the drive.
- Remove customer's materials.
- Place any DVD or MOD in the customer material box.
- Raise table to maximum height and return gantry to zero. (Check to verify tilt brackets fit.)

3-6 Power-down Console

NOTICE

Complete all work that requires system power "ON" before performing this process step. This may include: gantry base cover removal, raise or lowering the table, etc.

1. Click the [Service Desktop] icon.
2. Click [Shutdown].
3. When the Power Off message appears, wait two minutes, and then turn off console power switch.
4. Shut down gantry power switches.

Chapter 4- De-installation Power Down & LOTO

The mechanical de-installation team should not perform de-installation work on any CT systems, or around the PDU until implementing Lock-Out Tag-Out (LOTO) on the system.

4-1 A-1 Disconnect Lock Out Procedure

Some systems include a UPS system as an emergency power backup. If your system includes a UPS, review the power on/off section in the UPS service manual before starting any de-installation procedures.



WARNING

Lock Out – Tag Out is required before performing any de-installation work. If not completed you could be electrocuted causing death or severe injury.

1. Follow all GEH procedures for LOTO of the main disconnect.
2. Turn off PDU breakers.
3. Turn off UPS power and breakers.
4. Turn off gantry breakers.
5. Measure voltage at the PDU main terminal to confirm that A1 power is off.



NOTICE

Room warning lights leads could still carry voltage. Use caution when removing these wires at the PDU.

6. **Do not cut system cables.**

All items should be packed for reuse unless noted in writing from the PMI. Refer to packaging guidelines in this manual.

Chapter 5- Mechanical De-installation Overview

1. Locate the installation manual for the system being removed. If the paper manual is not found, locate the service CD and use the installation manual on the disk.
2. Review the installation procedures for each task prior to disassembly. The procedures in the de-installation manual are to be used as a guideline to disassembling the components. Detailed disassembly procedures are not required.
3. Gantry tilt brackets and gantry leveler shipping bolts are often located in the FE service cabinet. These items will be reused if the system is being reinstalled.
4. If a de-installation packaging kit was ordered, locate and unpack the items, which should include:
 - Dollies
 - Boxes
 - Packaging materials
 - Skids
 - Console skids
5. Prior to packing equipment in boxes, place the large boxes on skids and ensure the skid will fit through doorways.
6. A typical CT component de-installation can be packed into two large boxes. A third large box is included if necessary. There are also smaller boxes for packing service documentation.
7. **Do not cut system cables.** All cables should remain in tact unless the Project Manager of Installation has specified otherwise. Pack and return all installed system cables.
8. The console should be mounted onto the return-shipping skid, supplied. **Reminder: Do not lift the console onto the skid.** If this is a VCT, remove the hard drives and pack them in the 'supplies box' found in the console. CRT monitors should be packed into the large supplied boxes.
9. The table footswitch assembly is normally located in a bio-contaminated area. **Reminder: Wear and use appropriate BBP procedures when working in and around contaminated areas.**
10. Return the gantry and table floor leveling pads with the system. Use supplied cable ties to secure pads, and pack into plastic bag.
11. Use ties to secure the gantry top covers to the frame. Use the packaging supplies to protect the gantry side covers. Place the two side covers on top of a filled shipping box and strap to skid. Shock watch devices shipped with the de-installation kit can be installed on the boxes and all other returned equipment.
12. The HP table base covers should be wrapped with protecting material shipped with packaging supplies. Place the two covers onto of table and strap to secure.
13. If the system is being returned, pack all system pads, phantoms and accessories in a box and label as 'Accessories'.

Chapter 6- De-Install the Operator Console

6-1 Remove & Pack Monitors, SCIM Assembly and SCSI Tower

1. Power and signal cables stay with the console.
2. Pack monitors in supplied packaging.

6-2 Disassemble & Pack Keyboard Table Top Accessories

1. Pack keyboard assembly in supplied packaging.
2. Pack the mouse and trackball in a supplied box.
3. Pack the bar code reader if present in supplied box.

6-3 VCT ONLY- Remove Hard Drives

1. Remove the hard drives and place in the drives box located inside the console.
2. Replace the box with the drives back inside the console for shipment.
3. Replace all removed covers.

6-4 Remove Console Rear Cable Connections

- All LAN connections.
- All gantry signal cables
- Gantry fiber cable
- Console power cord
- Console #2 ground
- All other attached cables

6-5 Remove all Disconnected System Cables

6-6 Pack console for shipping

1. Pack console assembly in supplied packaging.
2. Remove the user desktop and arm assembly. Pack for shipment.
3. Load GO1 or later console onto supplied shipping skid, See Illustration 5-1.
Note: If loading a H1-H3 console use supplied dollies.
4. Pack all console accessory boxes on the console skid.



Illustration 5-1: GRE Console shown ready to ship.

All items should be packed for reuse unless noted in writing from the PMI. Refer to packaging guidelines in this manual.

Chapter 7- PDU and UPS De-installation

7-1 De-Install Power Distribution Unit

1. Remove PDU top and front covers
2. Remove all attached cables
 - a. Main feeders
 - b. Gantry and console power cables
 - c. System signal cables
 - d. System grounds
 - e. Warning light connection.



⚠ WARNING

Room warning lights leads could still carry voltage. Use caution when removing these wires at the PDU.

3. Remove PDU flexible cable conduit
4. If installed, remove all UPS connections
5. Confirm A1 power is disconnected. Remove remaining UPS connections
6. If installed, remove seismic mounting
7. Move PDU away from the wall
8. Install removed covers
9. Install supplied PDU packaging

7-2 De-Install Uninterrupted Power Supply (UPS)

1. Remove UPS front cover
2. Check that UPS is in the service mode
3. Turn off circuit breakers reinstall cover.



⚠ WARNING

Batteries are charged and voltage is present. Use caution when disconnecting all cables.

4. Disconnect all attached cables
5. Follow Powerware's instruction to disconnect the batteries for shipment.
6. Remove all flexible cable conduit connected to the UPS
7. If installed, remove seismic mounting brackets
8. Move UPS away from the wall

9. Install removed covers
10. Install supplied UPS packaging
11. Pack UPS cables and place with system.

All items should be packed for reuse unless noted in writing from the PMI. Refer to packaging guidelines in this manual.

Chapter 8- De-installation of the Table

1. Remove table footswitch cover. Save all hardware in supplied bags.
2. Remove table base covers. Save all hardware in supplied bags.
3. Remove table vertical side panels to access mounting anchors.
4. For GT tables, remove left side cradle cover. Install IMS transportation bracket.
5. Re-install left side cradle cover.
6. Remove table cables.
 - a. Power
 - b. Ground



⚠ WARNING

BIOHAZARDS COULD BE PRESENT WHEN REMOVING COMPONENTS AND COVERS, WHICH HAVE BEEN EXPOSED TO PATIENT CONTACT.

Follow your Organization's Blood Borne Pathogen procedures, and Policies. Follow all Personal Protective Equipment requirements.

7. Remove footswitch assembly
8. Reassembly footswitch and pack for shipping.
9. Using supplied hardware, attach table dolly
10. With the table dollies installed, un-tighten, and remove anchor nuts
11. Using the dolly, raise the table over the anchors and away from the anchors.
12. Collect the leveler pads and place in supplied bags. Place bag inside the table.
13. Reinstall all table side panels
14. Place protective covering on top of the cradle area, and around each end of the cradle.
15. Place the table base cover in supplied boxes. Place boxes on the cradle for shipment.
16. Collect any loose hardware and ship with system.
17. Cut all existing floor anchors unless noted in writing from the PMI.

All items should be packed for reuse unless noted in writing from the PMI. Refer to packaging guidelines in this manual.

Chapter 9- De-installation of the Gantry

1. Remove all 5 gantry covers.
 - a. Save all hardware in supplied bags.
 - b. Pack covers.

**⚠ WARNING**

BIOHAZARDS COULD BE PRESENT WHEN REMOVING COMPONENTS AND COVERS, WHICH HAVE BEEN EXPOSED TO PATIENT CONTACT.

Follow your Organization's Blood Borne Pathogen procedures, and Policies. Follow all Personal Protective Equipment requirements.

2. Remove all gantry base covers, and hardware.
3. If present, remove gantry rear cable cover assembly.
 - a. Save all hardware in supplied bags.
 - b. Pack covers.
4. Install tilt brackets on the gantry. The brackets are with supplied hardware.
5. Remove all gantry cables.
 - a. Power
 - b. Ground
 - c. System cables

⚠ CAUTION

Do not step-on or bend fiber optic cable. Bending or crushing will destroy cable.

6. Using hardware supplied with the de-installation kit, attach gantry dollies.

**NOTICE**

DO NOT PIN THE GANTRY FOR SHIPMENT. The gantry must have the ability to rotate during shipment or the impact from moving will damage components.

7. With the gantry dollies installed, un-tighten, and remove all anchor nuts
8. Using the dolly, raise the gantry over the anchors and away from the anchors.
 - a. Collect the leveler pads and place in supplied bags.
 - b. Place bag with gantry.
9. Reinstall gantry top covers and secure for shipping.
10. Place the gantry side cover in supplied boxes.

11. Collect any loose hardware and ship with system.
12. If present, pack the gantry rear cable cover assembly.
13. Cut all existing floor anchors, unless noted in writing from the PMI.

All items should be packed for reuse unless noted in writing from the PMI. Refer to packaging guidelines in this manual.

Chapter 10- System Options

10-1 Injectors

CAUTION

This procedure requires two people to safely complete this task. Review the installation instructions to familiarize your self with the assembly method. Use this knowledge in the disassembly of the component.

1. Disassemble the installed components starting at the injector head
2. Remove the arm from the ceiling mount or pedestal. **Use Caution, the arm will be spring-loaded.**
3. Use a ladder to access the ceiling mounted pedestal.
4. With one person supporting the pedestal, unbolt and lower to the pedestal to the floor.
5. Pull all cable and pack for shipment
6. Disconnect the power supply and the control panel
7. Appropriately pack all removed components.

10-1-1 Floor Mounted Injectors

1. Disconnect all attached cables
2. Remove the injector head from the stand
3. Disconnect the power supply and the control panel
4. Appropriately pack all removed components.

10-2 Monitors

CAUTION

This procedure requires two people to safely complete this task. Review the installation instructions to familiarize your self with the assembly method. Use this knowledge in the disassembly of the component.

1. Disassemble the installed components starting at the monitor.
2. Remove the arm from the ceiling mount. **Use Caution, the arm will be spring loaded.**
3. Use a ladder to access the ceiling mount.
4. With one person supporting the hardware, unbolt and lower to the floor.
5. Pull all cable and pack for shipment
6. Disconnect the power supply and the control panel
7. Appropriately pack all removed components.

10-3 RPM D4D - Respiratory Gating

Review the installation instructions to familiarize your self with the assembly method. Use this knowledge in the disassembly of the component.

1. Disassemble the installed components starting at the monitor assembly on the table.
2. Disconnect all installed connections and pack appropriately
3. Remove the flat tabletop and pack appropriately
4. Disconnect the ceiling cable and drape assembly.
5. Remove the external power supply assembly.
6. Remove the Varian host computer.
7. Remove all cables.
8. Appropriately pack all removed components.

All items should be packed for reuse unless noted in writing from the PMI. Refer to packaging guidelines in this manual.

Chapter 11- Cables, Accessories and Final Procedures

1. Remove all system cables that connect the Console, PDU, UPS, Table and Gantry sub-systems.
 - a. Coil cables so that they fit into supplied packaging.
 - b. Proper PPE required for this task.
 - c. Beware of biohazards when performing this task.
2. Collect all adjuster / lever pads and put in a bag. Store bag in accessories box.
3. Pack all remaining covers in supplied boxes.
4. Collect all service tools and pack in supplied boxes.
 - a. Phantoms
 - b. Phantom holder
 - c. Cover dollies front
 - d. Cover dollies rear
 - e. Tube hoist assembly
 - f. Service tools kit
 - g. Balance weights
5. Remove service documentation from service cabinet. Pack documentation in supplied boxes
6. Locate and pack head holders and foot extender
7. Remove FE service cabinet and pack for shipment
8. Pack the operators chair
9. Ensure all items are secure on pallets. Use Avis straps when necessary
10. Seal all boxes with tape.
11. Apply Shockwatch labels to appropriate areas. Follow the instructions enclosed with the labels.
12. Remove the LOTO lock. Apply the appropriate WARNING tag, which indicates the system has been disconnected and the PDU wires will be active if restored.

All items should be packed for reuse unless noted in writing from the PMI. Refer to packaging guidelines in this manual.

Appendix A

Returned Equipment Inventory List

Item #	Equipment / Material	How Packed for GS	How Packed for RR
1	CT Gantry	On Dollies	On Dollies
2	Table	On Dollies	On Dollies
3	PDU	Cardboard Box	Card board Box
4	Console w desktop	Shrink wrap on Shipping Skid	Shrink wrap on Shipping Skid
5	UPS	Shrink wrap on Skid	Shrink wrap on Skid
6	Service Cabinet	Shrink wrap on Skid	Shrink wrap on Skid
7	SCIM key board	Ship packed in box	Pack in group box
8	SCSI Tower	Ship packed in box	Pack in group box
9	Monitors (2)	Ship packed in box	Pack in group box
10	Phantom QA	Ship packed empty in box	Pack in group box
11	Phantom 48	Ship packed in box	Pack in group box
12	Phantom 30	Ship packed in box	Pack in group box
13	Mouse	Ship packed in box	Pack in group box
14	Trackball	Ship packed in box	Pack in group box
15	Bar Code Reader	Ship packed in box	Pack in group box
16	Chair	Return Optional	Return Optional
17	Head holder	Ship packed in box	Pack in group box
18	Head holder	Ship packed in box	Pack in group box
19	Foot extender	Ship packed in box	Pack in group box
20	Tech -Pubs	Ship packed in box	Pack in group box
21	System software	Ship packed in box	Pack in group box
22	All dollies	Ship packed in box	Pack in group box
23	Tube hoist asm	Ship packed in box	Pack in group box
24	Detector lift asm	Ship packed in box	Pack in group box
25	ECG option	Ship packed in box	Ship packed in box
26	D3D option	Ship packed in box	Ship packed in box
27	Flat table top	Ship packed in box	Pack in group box
28	Boom in room	Ship packed in box	Pack in group box
29	IVY Monitor	Ship packed in box	Ship packed in box
30	Cables	Ship packed in box	Pack in group box
31	All ceiling mounted equipment	Ship packed in box	Pack in group box
32	Foot switch assembly	Ship packed in box	Ship packed in box

Customer Items**Do Not Remove**

Item #	Equipment / Material
1	Customer pads & straps
2	Customer MOD's
3	Customer documents
4	Customer films
5	FE parts, tools, components
5	FE software

Goldseal Returned Packaging

Item #	Equipment / Material	Dolly / packaging descriptions	Differences
1	CT Gantry	Standard CT gantry dollies	Dollies are different for each gantry
2	Table	Standard CT table dollies	Dollies are different for each gantry
3	PDU	Cardboard Box	Standard PDU packaging
4	Console w desktop	Standard CT console Shipping Skid	Shrink wrap on Shipping Skid
5	UPS	Need a skid for this item	This item is heavy
6	Service Cabinet	Standard Skid	Shrink wrap on Skid
7	SCIM key board	Need one box with 4 dividers for desktop items	
8	SCSI Tower	See 7	
9	Monitors (2)	Need one box with 2 dividers	
10	Phantom QA	Ship packed empty in box	
11	Phantom 48	Pack wrapped in one box	
12	Phantom 30	Pack wrapped in one box	
13	Mouse	See 7	
14	Trackball	See 7	
15	Bar Code Reader	See 7	
16	Chair	Return Optional	Return Optional
17	Head holder	Pack wrapped in one box	
18	Head holder	Pack wrapped in one box	
19	Foot extender	Pack wrapped in one box	
20	Tech -Pubs	Ship packed in box	
21	System software LFC	Ship packed in box	
22	All dollies	Ship packed in box	Heavy
23	Tube hoist asm	Ship packed in box	
24	Detector lift asm	Ship packed in box	
25	ECG option	Ship packed in 2 boxes	
26	D3D option	Ship packed in 3 boxes	
27	Flat table top	Ship packed in box	Long item do not ding
28	Boom in room	Ship packed in 3 boxes	
29	IVY Monitor	Ship packed in 2 boxes	
30	Cables	Ship packed in 2 boxes	
31	All ceiling mounted equipment	Ship packed in 3 boxes	
32	Foot switch assembly	Ship packed in box	Bio Hazard

This page was intentionally left blank.

(c) 2008, General Electric Company.
GE Medical Systems, a General Electric Company, going to market as GE Healthcare.
3000 N. Grandview Boulevard
Waukesha, Wisconsin 53188
USA

www.gehealthcare.com



Imagination at work