

**ADMINISTRATIVE INFORMATION****FMI No: 22901****FMI Title: NGPDU SOFT-START RESISTOR REPLACEMENT**

FMI TYPE	CHARGE CODES	IMPLEMENTATION DATES
☐ Mandatory Safety	GEMS-AM: Code AM	Issue Date:
☒ Mandatory	GEMS-E: Code E	Close Date:
☐ Mandatory on Request	GEMS-A: Code A	

☒ Staged FMI	Customer List provided by FMI Administration. If you identify additional units within effectivity range, but they are not shown on the customer list, contact FMI Administration in your Pole.
☒	Parts and any special tools are furnished automatically. To address shortages and ordering exceptions, contact FMI Administration in your Pole.
☐	Parts Ordering: GEMS - Americas, contact your Region logistics coordinator GEMS - Europe, Order FMI kits from GPO/EPL Buc GEMS - Asia, Contact your region logistics coordinator
☐	No Parts are required to complete this FMI.
☐ Batch with	<u>FMI #</u> To minimize costs, combine implementation of this FMI with the FMIs listed for all customers where the effectivity is the same (Please see your customer lists).
☐ Non-Staged	No Customer List is provided. Each service location must identify affected units and place part orders as indicated in the Special Instructions.

Technical Support

GEMS - Americas	GEMS - Europe	GEMS - Asia
GE Medical Systems Insite Center (414) 524-5300 800-321-7973 Milwaukee, WI 53201, USA	GEMS - E / ESC BP34 F 78533 BUC CEDEX, FRANCE Tel: (33 1) 39 20 00 07 or ESC Fax: (33 1) 30 70 99 70	Asia Support Center Phone: 81-426-56-0033 Fax: 81-426-48-2905

Special Instructions:**Unit Tracked:****FMI Instructions Part No.: 2405311****Do NOT leave unsecured on site.**



Completion Certification Report

Mail To

GEMS - Americas

GE Medical Systems
Attn: FMI Admin. W-523
P.O. Box 414
Milwaukee, WI 53201

GEMS - Europe

GEMS - E / ESC
FMI Admin. 322
BP34
F 78533 Buc Cedex, France

Send to your
Country FMI or Service
Administrator

FMI No.: 22901 FMI Title: NGPDU SOFT-START RESISTOR REPLACEMENT

Customer Full Name: _____

Number and Street Address: _____

City and Country: _____

Country Code

Site Number

System I.D.

Unit Tracked:

Model Number	Serial Number	Compl. Date DD/MM/YY	FMI Comp. Code	Service Hours	
				Labor	Travel

FMI Completion Code:

- 1.) Modification Complete or previously modified.
- 2.) FMI not required per FMI instruction.
- 3.) Not modified. Unit in storage. FMI kit returned to stock.
- 4.) Modification refused or equipment altered. Customer signature and title required.
- 5.) Unit location unknown or unit not in area. Product Locator notified.

Customer Refusal or Equipment Altered - Modification Pending.

Installer
Comments _____

EFFECTIVITY:

The effective products are Hispeed QX/i, LightSpeed Pro16, LightSpeed 16 (H-Power 60), LightSpeed RT and LightSpeed Plus/Ultra which use NGPDU with the following serial numbers. (In the list, SN: Serial Number; Date: Manufactured Date; Serial numbers are listed in the ascending order, from top to bottom and from left to right.)

Note: When you spot the serial number of your site table, note the manufactured date.

<u>SN</u>	<u>Date</u>	<u>SN</u>	<u>Date</u>	<u>SN</u>	<u>Date</u>
50928HM1	09-Dec-02	56415HM3	26-Mar-03	57165HM3	22-Apr-03
50929HM3	09-Dec-02	56416HM1	26-Mar-03	57166HM1	22-Apr-03
50930HM5	09-Dec-02	56417HM9	26-Mar-03	57698HM3	22-Apr-03
50931HM7	09-Dec-02	56418HM7	26-Mar-03	57699HM1	22-Apr-03
52641HM8	09-Jan-03	56419HM5	23-Mar-03	57700HM7	22-Apr-03
52642HM6	09-Jan-03	56728HM9	01-Apr-03	57701HM5	22-Apr-03
53229HM1	17-Jan-03	56729HM7	01-Apr-03	57702HM3	22-Apr-03
53230HM9	17-Jan-03	56730HM5	01-Apr-03	57703HM1	22-Apr-03
53972HM6	12-Feb-03	56731HM3	01-Apr-03	57704HM9	22-Apr-03
53973HM4	19-Feb-03	56732HM1	08-Apr-03	57705HM6	22-Apr-03
54825HM5	03-Mar-03	56733HM9	08-Apr-03	57706HM4	22-Apr-03
54826HM3	03-Mar-03	56734HM7	08-Apr-03	57707HM2	22-Apr-03
54827HM1	03-Mar-03	56735HM4	08-Apr-03	57708HM0	22-Apr-03
54828HM9	03-Mar-03	56736HM5	31-Mar-03	57709HM8	22-Apr-03
55182HM0	04-Mar-03	56737HM0	04-Apr-03	57710HM6	22-Apr-03
55183HM8	04-Mar-03	56738HM8	04-Apr-03	58100HM9	28-Apr-03
55184HM6	04-Mar-03	56739HM6	04-Apr-03	58101HM7	28-Apr-03
55185HM3	04-Mar-03	56740HM4	04-Apr-03	58102HM5	28-Apr-03
55186HM1	10-Mar-03	56741HM2	07-Apr-03	58103HM3	28-Apr-03
55187HM9	10-Mar-03	56742HM0	07-Apr-03	58104HM1	28-Apr-03
55188HM7	10-Mar-03	56743HM8	07-Apr-03	58105HM8	28-Apr-03
55189HM5	10-Mar-03	56744HM6	07-Apr-03	58291HM6	07-May-03
55190HM3	07-Mar-03	56745HM3	07-Apr-03	58292HM4	07-May-03
55191HM1	07-Mar-03	57159HM6	08-Apr-03	58293HM2	07-May-03
56410HM4	26-Mar-03	57160HM4	08-Apr-03	58294HM0	07-May-03
56411HM2	26-Mar-03	57161HM2	08-Apr-03	58295HM7	07-May-03
56412HM0	26-Mar-03	57162HM0	08-Apr-03	58691HM7	09-May-03
56413HM8	26-Mar-03	57163HM8	18-Apr-03	58692HM5	09-May-03
56414HM6	26-Mar-03	57164HM6	22-Apr-03	58693HM3	09-May-03

<u>SN</u>	<u>Date</u>	<u>SN</u>	<u>Date</u>	<u>SN</u>	<u>Date</u>
58694HM1	10-May-03	59807HM8	07-Jun-03	61013HM9	17-Jun-03
58695HM8	10-May-03	59808HM6	07-Jun-03	61014HM7	17-Jun-03
59095HM0	18-May-03	59809HM4	07-Jun-03	61015HM4	18-Jun-03
59096HM8	18-May-03	60357HM1	10-Jun-03	61016HM2	18-Jun-03
59097HM6	18-May-03	60358HM9	10-Jun-03	61017HM0	18-Jun-03
59098HM4	18-May-03	60359HM7	10-Jun-03	61018HM8	18-Jun-03
59099HM2	18-May-03	60360HM5	10-Jun-03	61019HM6	18-Jun-03
59100HM8	10-Sep-03	60361HM3	10-Jun-03	61599HM7	24-Jun-03
59101HM6	10-Sep-03	60362HM1	10-Jun-03	61600HM3	24-Jun-03
59102HM4	18-May-03	60363HM9	11-Jun-03	61601HM1	24-Jun-03
59114HM9	18-May-03	60364HM7	11-Jun-03	61602HM9	24-Jun-03
59421HM8	26-May-03	60365HM4	11-Jun-03	61603HM7	24-Jun-03
59422HM6	26-May-03	60366HM2	11-Jun-03	61604HM5	24-Jun-03
59423HM4	26-May-03	60367HM0	11-Jun-03	61605HM2	24-Jun-03
59424HM2	26-May-03	60368HM8	11-Jun-03	61606HM0	08-Jul-03
59425HM9	26-May-03	60369HM6	11-Jun-03	61607HM8	24-Jun-03
59426HM7	26-May-03	60370HM4	11-Jun-03	61608HM6	08-Jul-03
59427HM5	26-May-03	60905HM7	14-Jun-03	61609HM4	08-Jul-03
59428HM3	26-May-03	60907HM3	14-Jun-03	61610HM2	08-Jul-03
59429HM1	26-May-03	60908HM1	14-Jun-03	61611HM0	08-Jul-03
59430HM9	26-May-03	60909HM9	14-Jun-03	61612HM8	08-Jul-03
59431HM7	26-May-03	60999HM0	23-Jun-03	61613HM6	08-Jul-03
59689HM0	26-May-03	61000HM6	23-Jun-03	61614HM4	08-Jul-03
59795HM5	02-Jun-03	61001HM4	23-Jun-03	61615HM1	08-Jul-03
59796HM3	02-Jun-03	61002HM2	23-Jun-03	61616HM9	08-Jul-03
59797HM1	02-Jun-03	61003HM0	14-Jun-03	61617HM7	08-Jul-03
59798HM9	02-Jun-03	61004HM8	14-Jun-03	61618HM5	08-Jul-03
59799HM7	02-Jun-03	61005HM5	16-Jun-03	61619HM3	08-Jul-03
59800HM3	03-Jun-03	61006HM3	14-Jun-03	61620HM1	08-Jul-03
59801HM1	03-Jun-03	61007HM1	14-Jun-03	61621HM9	08-Jul-03
59802HM9	04-Jun-03	61008HM9	16-Jun-03	61622HM7	10-Jul-03
59803HM7	04-Jun-03	61009HM7	16-Jun-03	61623HM5	10-Jul-03
59804HM5	04-Jun-03	61010HM5	17-Jun-03	61624HM3	10-Jul-03
59805HM2	07-Jun-03	61011HM3	17-Jun-03	62526HM9	09-Jul-03
59806HM0	07-Jun-03	61012HM1	17-Jun-03	62527HM7	09-Jul-03

<u>SN</u>	<u>Date</u>	<u>SN</u>	<u>Date</u>	<u>SN</u>	<u>Date</u>
62528HM5	11-Jul-03	63743HM9	05-Aug-03	64848HM5	19-Aug-03
62529HM3	11-Jul-03	63744HM7	05-Aug-03	64849HM3	19-Aug-03
62530HM1	11-Jul-03	63745HM4	06-Aug-03	64850HM1	19-Aug-03
62531HM9	11-Jul-03	63746HM2	06-Aug-03	64851HM9	19-Aug-03
62532HM7	14-Jul-03	63747HM0	06-Aug-03	64852HM7	19-Aug-03
62533HM5	14-Jul-03	63748HM8	06-Aug-03	64853HM5	19-Aug-03
62534HM3	14-Jul-03	64191HM0	05-Aug-03	64972HM3	20-Aug-03
62535HM0	14-Jul-03	64192HM8	08-Aug-03	64973HM1	20-Aug-03
62536HM8	14-Jul-03	64193HM6	08-Aug-03	64974HM9	20-Aug-03
62537HM6	14-Jul-03	64194HM4	08-Aug-03	64975HM6	20-Aug-03
63217HM4	22-Jul-03	64195HM1	08-Aug-03	65329HM5	25-Aug-03
63218HM2	22-Jul-03	64196HM9	08-Aug-03	65330HM3	25-Aug-03
63219HM0	22-Jul-03	64197HM7	12-Aug-03	65331HM1	25-Aug-03
63220HM8	22-Jul-03	64198HM5	12-Aug-03	65332HM9	30-Sep-03
63221HM6	22-Jul-03	64199HM3	12-Aug-03	65333HM7	23-Sep-03
63222HM4	22-Jul-03	64200HM9	12-Aug-03	65334HM5	23-Sep-03
63223HM2	22-Jul-03	64201HM7	12-Aug-03	65335HM2	25-Aug-03
63224HM0	22-Jul-03	64202HM5	12-Aug-03	65336HM0	25-Aug-03
63225HM7	22-Jul-03	64203HM3	12-Aug-03	65337HM8	25-Aug-03
63226HM5	22-Jul-03	64204HM1	12-Aug-03	65338HM6	26-Aug-03
63227HM3	23-Jul-03	64299HM1	10-Sep-03	65339HM4	26-Aug-03
63228HM1	23-Jul-03	64300HM7	23-Sep-03	65340HM2	26-Aug-03
63229HM9	23-Jul-03	64301HM5	23-Sep-03	65341HM0	26-Aug-03
63563HM1	29-Jul-03	64302HM3	13-Aug-03	65342HM8	27-Aug-03
63564HM9	29-Jul-03	64303HM1	13-Aug-03	65343HM6	27-Aug-03
63665HM4	29-Jul-03	64304HM9	13-Aug-03	65344HM4	27-Aug-03
63666HM2	29-Jul-03	64305HM6	15-Aug-03	65729HM6	24-Oct-03
63667HM0	29-Jul-03	64306HM4	15-Aug-03	65730HM4	24-Oct-03
63668HM8	29-Jul-03	64307HM2	15-Aug-03	65731HM2	24-Oct-03
63669HM6	29-Jul-03	64308HM0	15-Aug-03	65732HM0	03-Sep-03
63670HM4	29-Jul-03	64309HM8	20-Aug-03	65733HM8	03-Sep-03
63671HM2	29-Jul-03	64310HM6	20-Aug-03	65734HM6	03-Sep-03
63672HM0	29-Jul-03	64311HM4	10-Sep-03	65735HM3	03-Sep-03
63673HM8	29-Jul-03	64312HM2	10-Sep-03	65736HM1	04-Sep-03
63742HM1	05-Aug-03	64847HM7	19-Aug-03	65737HM9	04-Sep-03

<u>SN</u>	<u>Date</u>	<u>SN</u>	<u>Date</u>	<u>SN</u>	<u>Date</u>
65738HM7	04-Sep-03	67952HM2	26-Sep-03	69893HM6	21-Oct-03
65903HM7	05-Sep-03	67953HM0	26-Sep-03	69894HM4	21-Oct-03
65904HM5	05-Sep-03	67954HM8	26-Sep-03	69895HM1	21-Oct-03
65905HM2	05-Sep-03	67955HM5	30-Sep-03	69896HM9	24-Oct-03
65906HM0	05-Sep-03	67956HM3	25-Sep-03	69897HM7	21-Oct-03
66280HM9	11-Sep-03	68058HM7	23-Sep-03	69898HM5	24-Oct-03
66281HM7	11-Sep-03	68952HM1	21-Oct-03	69899HM3	27-Oct-03
66282HM5	11-Sep-03	68953HM9	21-Oct-03	69900HM9	28-Oct-03
66283HM3	12-Sep-03	68954HM7	21-Oct-03	70542HM6	27-Oct-03
66284HM1	12-Sep-03	68955HM4	21-Oct-03	70543HM4	28-Oct-03
66285HM8	12-Sep-03	68956HM2	21-Oct-03	70544HM2	28-Oct-03
66286HM6	12-Sep-03	68957HM0	21-Oct-03	70545HM9	28-Oct-03
66980HM4	18-Sep-03	68958HM8	30-Sep-03	70546HM7	29-Oct-03
66981HM2	18-Sep-03	68959HM6	30-Sep-03	70547HM5	29-Oct-03
66982HM0	18-Sep-03	68960HM4	30-Sep-03	70548HM3	29-Oct-03
66983HM8	18-Sep-03	68961HM2	30-Sep-03	70559HM0	28-Oct-03
66984HM6	22-Sep-03	69233HM5	14-Oct-03	70560HM8	27-Oct-03
66985HM3	22-Sep-03	69234HM3	14-Oct-03	70561HM6	27-Oct-03
66986HM1	22-Sep-03	69235HM0	14-Oct-03	71056HM6	02-Nov-03
66987HM9	22-Sep-03	69236HM8	14-Oct-03	72021HM9	11-Nov-03
67443HM2	16-Sep-03	69237HM6	14-Oct-03	73398HM0	25-Nov-03
67445HM7	16-Sep-03	69238HM4	14-Oct-03	73399HM8	25-Nov-03
67446HM5	17-Sep-03	69239HM2	14-Oct-03	73400HM4	25-Nov-03
67447HM3	17-Sep-03	69240HM0	15-Oct-03	73401HM2	25-Nov-03
67448HM1	23-Sep-03	69241HM8	15-Oct-03	73999HM5	30-Nov-03
67449HM9	23-Sep-03	69242HM6	14-Oct-03	74000HM1	30-Nov-03
67450HM7	23-Sep-03	69243HM4	14-Oct-03	75028HM1	07-Dec-03
67451HM5	23-Sep-03	69244HM2	15-Oct-03	75029HM9	07-Dec-03
67452HM3	17-Sep-03	69245HM9	15-Oct-03	75030HM7	07-Dec-03
67453HM1	22-Sep-03	69246HM7	14-Oct-03	75031HM5	07-Dec-03
67947HM2	23-Sep-03	69247HM5	15-Oct-03	75032HM3	07-Dec-03
67948HM0	25-Sep-03	69889HM4	28-Oct-03	75033HM1	07-Dec-03
67949HM8	26-Sep-03	69890HM2	28-Oct-03	75288HM1	09-Dec-03
67950HM6	26-Sep-03	69891HM0	28-Oct-03	75289HM9	09-Dec-03
67951HM4	25-Sep-03	69892HM8	21-Oct-03	75290HM7	09-Dec-03

<u>SN</u>	<u>Date</u>	<u>SN</u>	<u>Date</u>	<u>SN</u>	<u>Date</u>
75626HM2	14-Dec-03	75630HM4	14-Dec-03		

PURPOSE:

Solve the issue that the old soft-start resistors frequently burn and improve the reliability of the NGPDU. Refer to PSR13012939/PQR13016041

RELATED FMI'S:

None

MANPOWER REQUIREMENTS:

One field engineer, one hour labor plus travel

SPECIAL TOOLS & TEST EQUIPMENT

None

FURNISHED MATERIALS

This FMI kit (2405313) consists of the following:

ITEM	PART NO.	DESCRIPTION	QTY.	COMMENT
1	2391041	New Soft-start Resistor	1	
2	2405311	FMI22901 Instruction	1	this manual

PROCEDURE



WARNING For systems powered from a Power Distribution Box which has a lockout capability, Turn OFF and lock out the system power at the box during the replacement procedure.



CAUTION ELECTRICAL HAZARD!
Please be assured that you have understood all the safety guidelines that you have to follow during your work.

Please Turn OFF the Main Power Distribution Supply to CT system and Switch OFF the CB1 and CB5 during your FMI work.

1 PREPARATION

Understand and follow the Safety Guidelines Manual.

- 1.) Switch OFF the OC main power.
- 2.) Switch OFF the breaker in the Power Distribution Box (hospital supply).

Note: If there is UPS, Switch OFF the UPS.

- 3.) Remove the Top Cover.
- 4.) Remove the Rear Cover.
- 5.) Remove the Plastic cover.
- 6.) Switch OFF the CB1 and CB5 breakers in the NGPDU.

2 RESISTOR REPLACEMENT


- 1.) Disassembly the old Resistor Assy (P/N: 2351496).
 - a.) Mark and disconnect the three plastic connectors.
 - b.) Unscrew the six screws.
 - c.) Remove the whole old Resistor Assy.

Figure 1 - Old Resistor Assy



- 2.) Install the new Resistor Assy (P/N: 2391041) to the NGPDU.
 - a.) Fix the Resistor Assy with the six old screws(a).
In order to facilitate the installation, loosen the screws(b) at the two ends of the three green resistors.
 - b.) Connect the three plastic connectors with their counterparts.

Figure 2 - New Resistor Assy

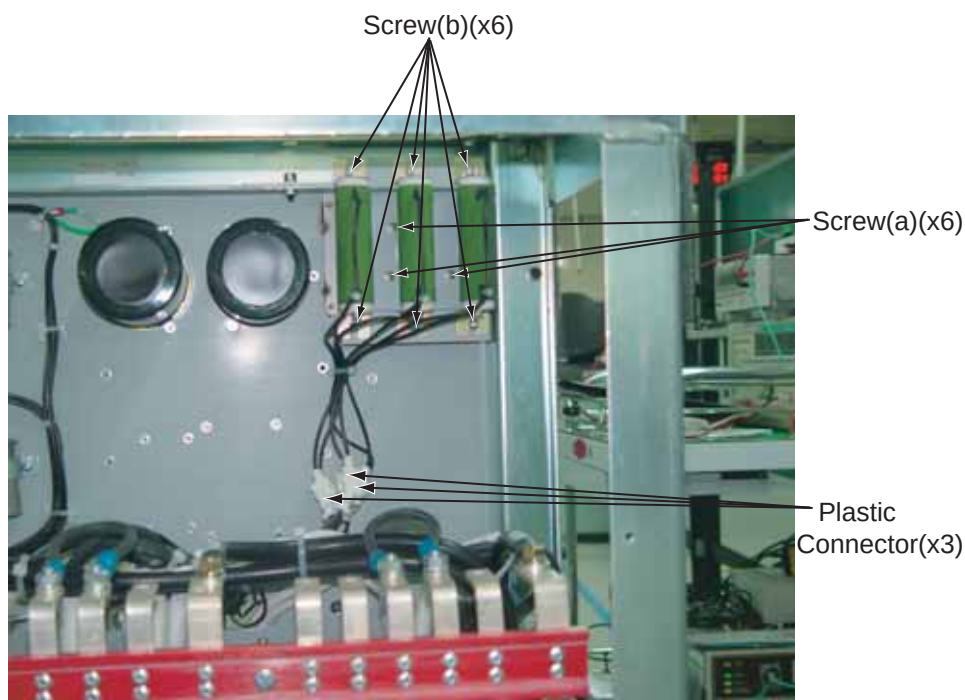

- 3.) Take a complete inspection to the Mechanical and Electrical connection.

3 FINAL CHECK

- 1.) Switch ON the Power Distribution Box.
- 2.) Switch ON the CB1 and CB5 breakers in the NGPDU.
- 3.) Switch ON the main switch of the OC. Verify that the system starts up with no error.
- 4.) Restart the system. Perform several scan with the system's maximal KV or mA.
Verify that the system can operate smoothly.

4 COMPLETING THE FMI

- 1.) Dispose of all removed / disused materials in accordance with your local policies for scrap material.
- 2.) Install the NGPDU covers.
- 3.) For countries equipped with Field MUST, GE CARES, or similar on-line systems; use your laptop computer to file your service report.
- 4.) For other countries; fill in the FMI completion certificate report and send them to your FMI coordinator.

