

**ADMINISTRATIVE INFORMATION****FMI No: 25313****FMI Title: LS2002 v305 LightSpeed Ultra/Plus MDAS Software FMI**

FMI TYPE	CHARGE CODES	IMPLEMENTATION DATES
☒ Mandatory Safety ☐ Mandatory ☐ Mandatory on Request	GEMS-AM: GEMS-E: GEMS-A:	Issue Date: Close Date:

☐ Staged FMI	Customer List provided by FMI Administration. If you identify additional units within effectivity range, but they are not shown on the customer list, contact FMI Administration in your Pole. ☐ Parts and any special tools are furnished automatically. To address shortages and ordering exceptions, contact FMI Administration in your Pole. ☐ Parts Ordering: GEMS - Americas, contact your Region logistics coordinator GEMS - Europe, Order FMI kits from GPO/EPL Buc GEMS - Asia, Contact your region logistics coordinator ☐ No Parts are required to complete this FMI.
☐ Batch with	To minimize costs, combine implementation of this FMI with the FMIs listed for all customers where the effectivity is the same (Please see your customer lists).
☐ Non-Staged	No Customer List is provided. Each service location must identify affected units and place part orders as indicated in the Special Instructions.

Technical Support

GEMS - Americas	GEMS - Europe	GEMS - Asia
GE Medical Systems Insite Center (414) 524-5300 800-321-7973 Milwaukee, WI 53201, USA	GEMS - E / ESC BP34 F 78533 BUC CEDEX, FRANCE Tel: (33 1) 39 20 00 07 or ESC Fax: (33 1) 30 70 99 70	Asia Support Center Phone: 81-426-56-0033 Fax: 81-426-48-2905

Special Instructions:**Unit Tracked:****FMI Instructions Part No.: 2358471-100****Do NOT leave unsecured on site.**



Completion Certification Report

Mail To

GEMS - Americas

GE Medical Systems
Attn: FMI Admin. W-523
P.O. Box 414
Milwaukee, WI 53201

GEMS - Europe

GEMS - E / ESC
FMI Admin. 322
BP34
F 78533 Buc Cedex, France

Send to your
Country FMI or Service
Administrator

FMI No.: 25313 FMI Title: LS2002 v305 LightSpeed Ultra/Plus MDAS Software FMI

Customer Full Name: _____

Number and Street Address: _____

City and Country: _____

Country Code

Site Number

System I.D.

Unit Tracked:

Model Number	Serial Number	Compl. Date DD/MM/YY	FMI Comp. Code	Service Hours	
				Labor	Travel

FMI Completion Code:

- 1.) Modification Complete or previously modified.
- 2.) FMI not required per FMI instruction.
- 3.) Not modified. Unit in storage. FMI kit returned to stock.
- 4.) Modification refused or equipment altered. Customer signature and title required.
- 5.) Unit location unknown or unit not in area. Product Locator notified.

**Customer Refusal or Equipment Altered
- Modification Pending.**

Installer
Comments

Effectivity

LightSpeed Ultra and Plus (MDAS) systems with the following Gantry model numbers are affected:

- 2281177
- 2281177-2
- 2281177-3
- 2281177-4

To more effectively track the location of the systems for this FMI, use any of the gantry models listed above, having a mfg date of December 6, 2002, or earlier.



NOTICE Do NOT load this software on Discovery LS Systems.

Purpose

The purpose of this FMI is to fix several software bugs and to introduce new applications to the customer, such as Cardiac Burst, Perfusion2, VariViewer, SmartStep, Auto mA and 440 mA Scan option.

This FMI will also fix FPR 10210992 (Mechanical Characterization) and FPR 1028447 (Biopsy screen displays wrong table location).

Furnished Materials

FMI 25313 Kit (Part Number 2361709) includes the following items:

Qty	Description	Part #
1	FMI Document	2358471-100
1	LS2002 v305 SW Set (MDAS) (see details, below)	2361705

The LS2002 SW Set (MDAS) (P/N 2361705) includes the following items:

Qty	Description	Part #
1	LS2002 v305 H3 M4 OS and Apps (MDAS) (see details, top of next page)	2360606-2
1	LS2002 Software Installation Manual	2348269-100
1	Application Tips and Workarounds	2296447-100
1	3.0 LightSpeed Plus Tech Ref Manual	2302180-100
1	CT LightSpeed Ultra CD-ROM	2304119-247
1	Serial #, Model #	46-302200P2
1	Options MOD (see details, on next page)	2361704

The LS2002 MDAS M4 Software Kit (P/N 2360606-2) includes the following items:

Qty	Description	Part #
1	Operating System Disk 1	2356015
1	Operating System Disk 2	2356016
1	Operating System Disk 3	2356017
1	Applications Software v305, Disk 1	2360608-2
1	Applications Software v305, Disk 2	2360609-2
1	LightSpeed Ultra/Plus LS2002 M4 Patch C	2366813

The Options MOD (P/N 2361704) includes the following:

Qty	Description	Part #
1	440 mA Scan Option	2345411
1	Auto mA Option	2356009
1	VariViewer Option*	2347142

* Customer must have 1.5GB of memory to load VariViewer Option. To confirm system has 1.5GB of memory, do the following:

1.) Open a shell

2.) Type: `hinv` ENTER

Look for main memory size of: 1536Mbytes

Required Tools

Qty	Description	Part #
1	LightSpeed Ultra/Plus (MDAS) Class C	2360598-2
1	LightSpeed Ultra/Plus (MDAS) Class M	2360596-2

Note: The `showprods` for Class C and Class M are as follow:


- Class C: `1100.x_CLS2002M4`
- Class M: `1200.x_MLS2002M4`

Note that `showprods` does not match a specific software version, because the Class C and Class M software for LightSpeed Plus (S/W v2.3) are the same as the Class C and Class M software for LightSpeed Ultra (S/W v305).

People Power

5 hours labor and 1 hour of travel.

Steering Guide

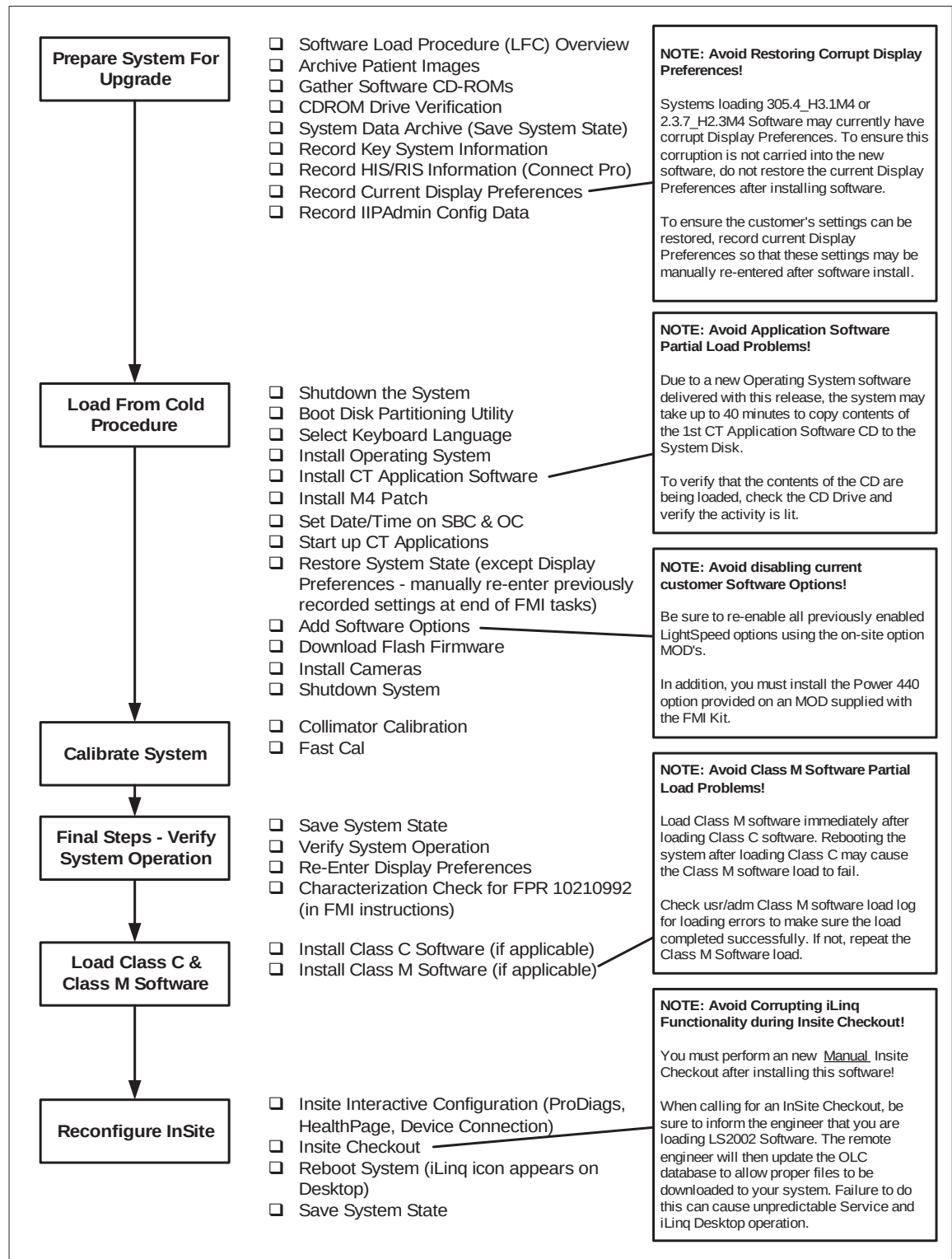

Procedure

1 Archive Patient Images

Verify that the customer has archived all patient data.

2 Software Load Procedures (LFC)

Perform the Load From Cold (LFC) procedure as detailed in Chapter 1 of *LS2002 Software Installation Procedures* (2348269-100).



NOTICE It is required that a “Restore STATE” be performed following the load of Class C and Class M software. To facilitate this, it is recommended that the Class C and Class M be loaded following the Applications Load (Section 3.6).

3 Calibration

- 1.) Perform Collimator Cal, by selecting UTILITIES on the bottom of the left monitor, then selecting COLLIMATOR CAL from the displayed selections.
- 2.) Perform FASTCAL by selecting DAILY PREP on the bottom of the left monitor, then selecting FASTCAL from the displayed selections.



NOTICE You may get a pop-up indicating: “System detects an alignment or change was performed on the System”. If this happens, press the Ok button and continue with the FASTCAL.

4 System Checkout and Save State

- 1.) Refer to the Software Installation Procedures Manual, and perform a Save System State. This will save the new calibrations as well as the customer and service protocols.
- 2.) Refer to the Software Installation Procedures Manual, and perform the system performance checks as described in the manual. Ensure all networking and filming processes are still functional.
- 3.) Complete and return Red Installation Product Locator Cards.

5 Characterization Check

5.1 Introduction

Use of reconfig in the prior software version, regardless of database reset selections, will delete current mechanical characterization files for tilt, elevation, and cradle longitudinal. These files' contents are replaced with default values.

The issue has been resolved in LS2002 software; however, the original system characterization files may have been lost and replaced with default values.

5.2 Procedure

- 1.) If not already completed, restore System State from the State MOD.
- 2.) Open a UNIX shell from the Tool Chest.
- 3.) At the prompt, enter the following command:

```
diff /usr/g/fw/char/tilt.def /usr/g/fw/char/tilt.char
```
- 4.) If any number is returned prior to receiving the command prompt, the two files are different. The implication is that this system is running the proper characterization files and no further action is required. Skip to step 7.
- 5.) If the prompt returns and no numbers are displayed prior to receiving the prompt, no difference exists between the files. The default characterization files are being used by the system. Proceed to step 6.
- 6.) Referring to the *System Service Manual* as needed, perform Elevation, Cradle and Tilt Mechanical Characterizations.
- 7.) Upon completion of step 6, save System State to MOD.

This Page Intentionally Blank.