

**ADMINISTRATIVE INFORMATION****FMI No: 25315****FMI Title: LightSpeed 16 M4 Software Upgrade**

FMI TYPE	CHARGE CODES	IMPLEMENTATION DATES
☒ Mandatory Safety	GEMS - AM:	Issue Date:
☐ Mandatory	GEMS - E:	
☐ Mandatory on Request	GEMS - A:	Close Date:

☒ Staged FMI	Customer List provided by FMI Administration. If you identify additional units within effectivity range, but they are not shown on the customer list, contact FMI Administration in your Pole.
☒	Parts and any special tools are furnished automatically. To address shortages and ordering exceptions, contact FMI Administration in your Pole.
☐	Parts Ordering: GEMS - Americas, Contact your Region logistics coordinator GEMS - Europe, Order FMI kits from GPO/EPL Buc GEMS - Asia, Contact your region logistics coordinator
☐	No Parts are required to complete this FMI.
☐ Batch with	_____
	To minimize costs, combine implementation of this FMI with the FMI's listed for all customers where the effectivity is the same (Please see your customer lists).
☐ Non-Staged	No Customer List is provided. Each service location must identify affected units, and place part orders as indicated in the Special Instructions.

Technical Support

GEMS - Americas	GEMS - Europe	GEMS - Asia
GE Medical Systems Insite Center (414) 524-5300 800-321-7973 Milwaukee, WI 53201, USA	GEMS - E / ESC BP34 F 78533 BUC CEDEX, FRANCE Tel: (33 1) 39 20 00 07 or ESC Fax: (33 1) 30 70 99 70	Asia Support Center Phone: 81-426-56-0033 Fax: 81-426-48-2905

Special Instructions:**Unit Tracked:****FMI Instructions Part No.: 2369739-100****Do NOT leave unsecured on site**

FMI Completion Certification Report

Mail To

GEMS - Americas

GEMS - Europe

G.E. Medical Systems
Attn: FMI Admin. W-523
P.O. Box 414
Milwaukee, WI 53201

GEMS - Europe
FMI Admin. 322
BP34
F78533 Buc Cedex, France

Send to your
Country FMI or Service
Administrator

FMI No.: 25315

FMI Title: LightSpeed 16 M4 Software Upgrade

Customer Full Name:

Number and Street Address:

City and Country:

Country Code

Site Number

System I.D.

Unit Tracked:

Model Number	Serial Number	Compl. Date DD/MM/YY	FMI Comp. Code	Service Hours	
				Labor	Travel

FMI Completion Code

1. Modification Complete or previously modified.
2. FMI not required per FMI instruction.
3. Not modified. Unit in storage. FMI kit returned to stock.
4. Modification refused or equipment altered. Customer signature and title required.
5. Unit location unknown or unit not in area. Product Locator notified.

Customer Refusal or Equipment Altered - Modification Pending.

Installer
Comments

Effectivity

The following CT Equipment is affected:

LightSpeed 16 systems with Gantry number 2339985 shipped prior to December 31, 2002

Purpose

The purpose of this FMI is to upgrade LightSpeed 16 CT scanners to the latest version software. This software includes quality fixes, and 3 new features (Direct3D, Variviewer, and Large Image Series). It also addressed FPR 10210992 (Mechanical Characterization) and FPR 1028447 (Biopsy Screen displays wrong table location).

Furnished Materials

Qty	Part number	Description
1	2369739-100	FMI document
1	2364535	LightSpeed 16 Apps Software Collector containing:
1	2364553	Application Disk #1
1	2364566	Application Disk #2
1	2356012	LightSpeed 16 OS Software Collector containing:
1	2356015	OS Disk #1
1	2356016	OS Disk #2
1	2325017	OS Disk #3
1	2348766-247	CT learning / Ref Guide
1	2348269-100	LS2002 Software Install Doc
	2352974-100	Apps Tips and Workarounds
	One Option MOD containing	
1	2347142	Variviewer Software Option
1	2268937	Direct 3D Software Option
1	2362597	Large Image Series Option


Tools Required

Qty	Description	Part #
1	LightSpeed 16 Class C Software	2364524
1	LightSpeed 16 Class M Software	2364525

People Power

5 hours labor and 1hours of travel

Steering Guide

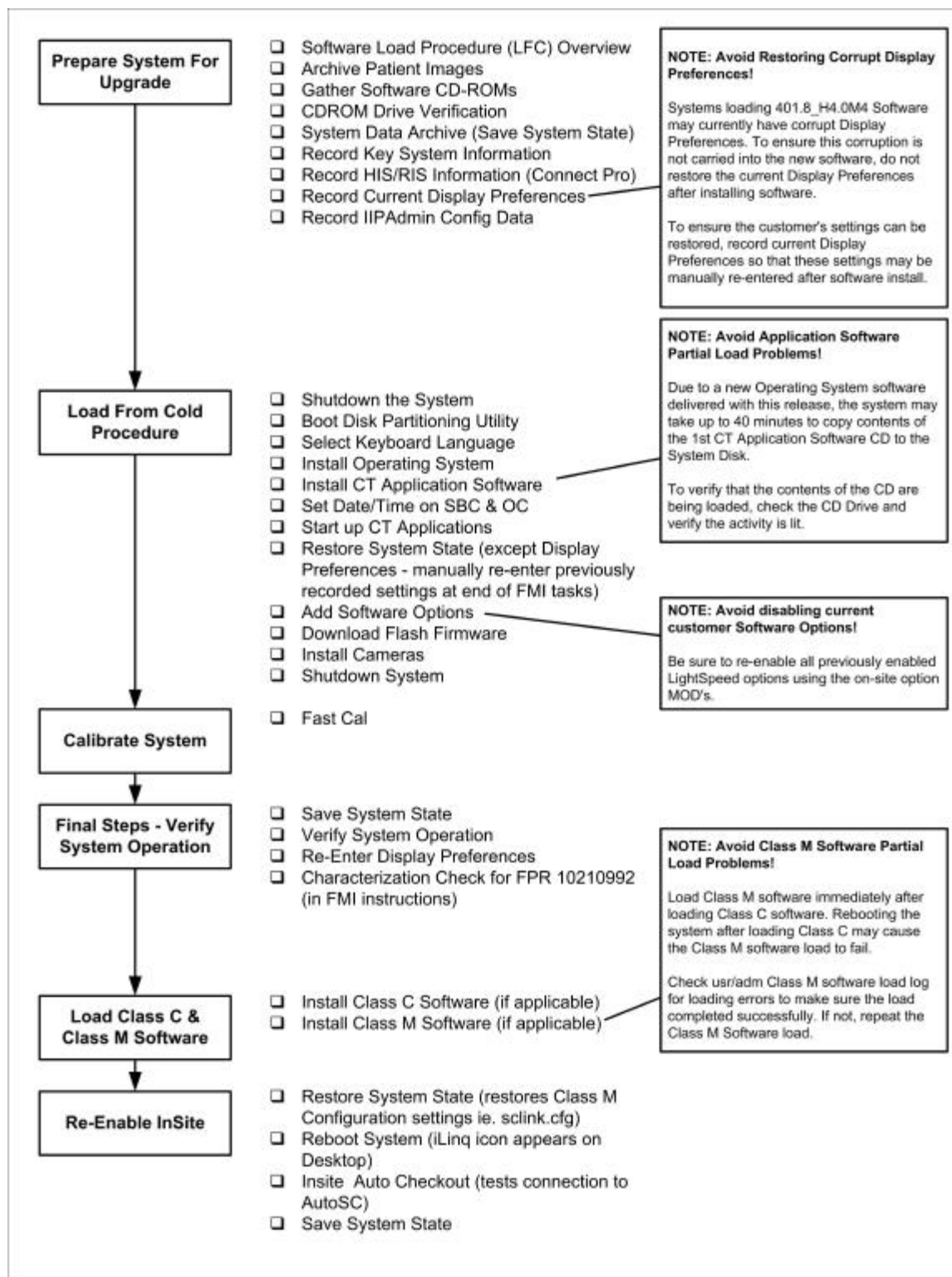

Procedure

1 Archive Patient Images

Verify that the customer has archived all patient data.

2 Software Load Procedures (LFC)

Perform the Load From Cold (LFC) procedure as detailed in Chapter 1 of *LS2002 Software Installation Procedures* (2348269-100).



NOTICE:

It is required that a “Restore STATE” be performed following the load of Class C and Class M software. To facilitate this, it is recommended that the Class C and Class M be loaded following the Applications Load.

3 Calibration

1. Perform Collimator Cal, by selecting UTILITIES on the bottom of the left monitor, then selecting COLLIMATOR CAL from the displayed selections.
2. Perform FASTCAL by selecting DAILY PREP on the bottom of the left monitor, then selecting FASTCAL from the displayed selections.



NOTICE:

You may get a pop-up indicating: “System detects an alignment or change was performed on the System”. If this happens, press the Ok button and continue with the FASTCAL.

4 System Checkout and Save State

- 1.) Refer to the Software Installation Procedures Manual, and perform a Save System State. This will save the new calibrations as well as the customer and service protocols.
- 2.) Refer to the Software Installation Procedures Manual, and perform the system performance checks as described in the manual. Ensure all networking and filming processes are still functional.
- 3.) Complete and return Red Installation Product Locator Cards.

5 Characterization Check

5.1 Introduction

Use of reconfig in the prior software version, regardless of database reset selections, will delete current mechanical characterization files for tilt, elevation, and cradle longitudinal. These files' contents are replaced with default values.

The issue has been resolved in LS2002 software; however, the original system characterization files may have been lost and replaced with default values.

5.2 Procedure

- 1.) If not already completed, restore System State from the State MOD.
- 2.) Open a UNIX shell from the Tool Chest.
- 3.) At the prompt, enter the following command:

```
diff /usr/g/fw/char/tilt.def /usr/g/fw/char/tilt.char
```
- 4.) If any number is returned prior to receiving the command prompt, the two files are different. The implication is that this system is running the proper characterization files and no further action is required. Skip to step 7.
- 5.) If the prompt returns and no numbers are displayed prior to receiving the prompt, no difference exists between the files. The default characterization files are being used by the system. Proceed to step 6.
- 6.) Referring to the System Service Manual as needed, perform Elevation, Cradle and Tilt Mechanical Characterizations.
- 7.) Upon completion of step 6, save System State to MOD.

Finishing Up

Complete and submit FMI closure forms as required by local practice.