

**ADMINISTRATIVE INFORMATION****FMI No: 25401A****FMI Title: VCT Cardiac Neuro and Auto Filter Options Installation**

FMI TYPE	CHARGE CODES	IMPLEMENTATION DATES
☐ Mandatory Safety	GEMS - AM:	Issue Date:
☒ Mandatory	GEMS - E:	
☐ Mandatory on Request	GEMS - A:	Close Date:

☒ Staged FMI	Customer List provided by FMI Administration. If you identify additional units within effectivity range, but they are not shown on the customer list, contact FMI Administration in your Pole.
☒	Parts and any special tools are furnished automatically. To address shortages and ordering exceptions, contact FMI Administration in your Pole.
☐	Parts Ordering: GEMS - Americas, Contact your Region logistics coordinator GEMS - Europe, Order FMI kits from GPO/EPL Buc GEMS - Asia, Contact your region logistics coordinator
☐	No Parts are required to complete this FMI.
☒ Batch with	FMI 25394 To minimize costs, combine implementation of this FMI with the FMI's listed for all customers where the effectivity is the same (Please see your customer lists).
☐ Non-Staged	No Customer List is provided. Each service location must identify affected units, and place part orders as indicated in the Special Instructions.

Technical Support

GEMS - Americas	GEMS - Europe	GEMS - Asia
GE Medical Systems Online Center 800-321-7937 Milwaukee, WI 53201, USA	GEMS - E / ESC BP34 F 78533 BUC CEDEX, FRANCE Tel: (33 1) 39 20 00 07 or ESC Fax: (33 1) 30 70 99 70	Asia Support Center Phone: 81-426-56-0033 Fax: 81-426-48-2905

Special Instructions:**This FMI applies ONLY to the following countries:**

- **Americas:** US and Latin America only
- **Europe:** All countries
- **Asia:** All countries excluding Korea, and China

If customer site is in Canada, Korea or China, do not load this FMI.**Unit Tracked:****FMI Instructions Part No.: 5224575-1EN****Do NOT leave unsecured on site**

FMI Completion Certification Report

Mail To

GEMS - Americas	GEMS - Europe	
GE Medical Systems Attn: FMI Admin. W-523 P.O. Box 414 Milwaukee, WI 53201, USA	GEMS - E / ESC BP34 F 78533 BUC CEDEX, FRANCE Tel: (33 1) 39 20 00 07 or ESC Fax: (33 1) 30 70 99 70	Asia Support Center Phone: 81-426-56-0033 Fax: 81-426-48-2905

FMI No.: 25401A

FMI Title: VCT Cardiac Neuro and Auto Filter Options Installation

Customer Full Name:

Number and Street Address:

City and Country:

Country Code

Site Number

System I.D.

Unit Tracked:

Model Number	Serial Number	Compl. Date DD/MM/YY	FMI Comp. Code	Service Hours	
				Labor	Travel

FMI Completion Code

1. Modification Complete or previously modified.
2. FMI not required per FMI instruction.
3. Not modified. Unit in storage. FMI kit returned to stock.
4. Modification refused or equipment altered. Customer signature and title required.
5. Unit location unknown or unit not in area. Product Locator notified.

Customer Refusal or Equipment Altered - Modification Pending.

Installer
Comments

Revision History

Revision	Date	Reason for Change
1	May 21, 2007	Initial Release

Effectivity

The FMI effectivity includes CT VCT64, VCT32 and VCTSTD systems. The following are tracked gantry model numbers.

Gantry Model Numbers:

- 5124069
- 5124069-2
- 5124069-5
- 5124069-6
- 5126093

This FMI is **ONLY** for VCT Systems in the following countries:

Americas: **Americas: US and Latin America only**

Europe: **All countries in Europe**

Asia: **All countries excluding Korea, and China**

Purpose

The purpose of this FMI is to install the following options made available in VCT Cardiac 07MW18.x software release (FMI 25394):

- Neuro Filters Software Option
- Auto Filter and Transfer Software Option

Prerequisite FMIs

FMI 25394.

Furnished Materials

FMI 25401A, Kit Part Number 5224575		
Qty	Description	Part #
1	FMI 25401A Document	5224575-1EN
1	Neuro Filters SW Option and Auto Filter and Transfer SW Option collector	5212526

Personnel Requirements

Travel Time: 1 hr

Labor Time: 0.5 hr

Procedure

1. Pre-requisite:

VCT Cardiac Software 07MW18.x (FMI 25393)

To check the software version:

- ☐ Open unix shell
 - o As a CT user
- ☐ Type **swinfo** <Enter>

2. Installation Options

- This section provides step-by-step detailed instructions on loading the options. If you are familiar with loading options from disk, and do not require detailed instructions, load the options now and continue with Section 3, 'Verification of the Options Install and final checks'.
- Insert Neuro Filter Option DVD-RAM (Part # 5199402) in the SCSI Tower DVD RAM drive.
- With Applications up, select the [**Service**] icon.
- Select [**Configuration**].
- Select [**Install Options**]. A screen similar to the one below appears (the screen will look different based on the options installed):

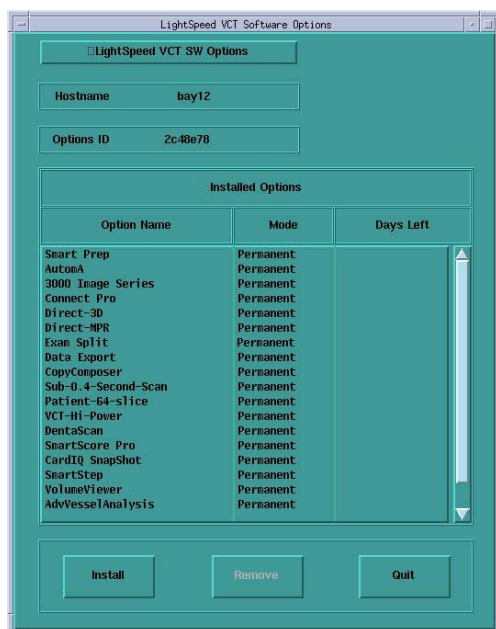


Figure 1 – Install Options Screen

6. Select [**Install**] . The Select Mechanism box appears.

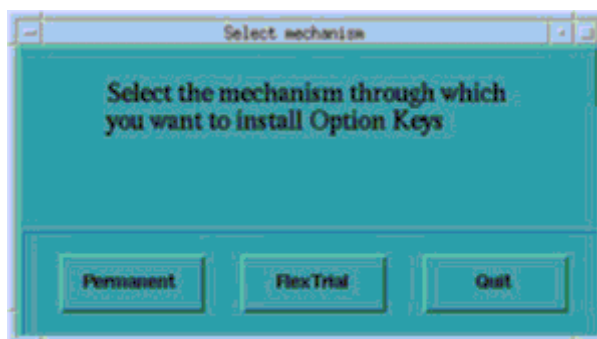


Figure 2 - Select Mechanism Box

7. Select [**Permanent**] . The Select Device box appears.



Figure 3 - Select Device Box

8. Select [**MEDIA**] . The Select Device box appears.

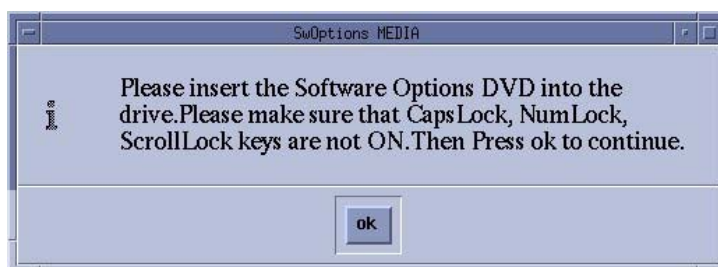


Figure 4 - Insert Option DVD Screen

9. Select [**ok**] for DVD insertion. The install options screen appears with the new installed option added to the list of options

10. Option selected will be shown in the screen

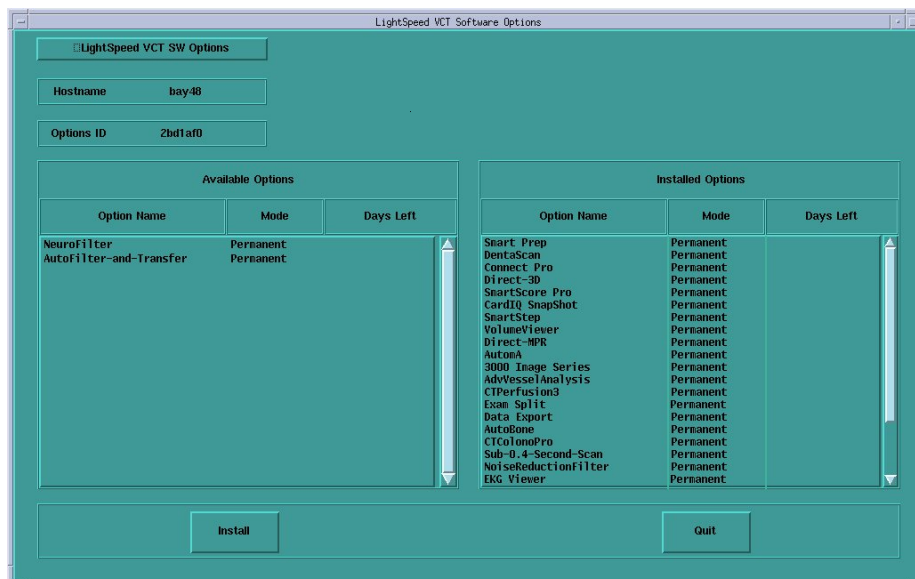


Figure 5 - Options Selected Example

11. Select both the Neuro Filter and auto Filter & Transfer options and select [Install]

The following screens appears one after the other indicating that options are installed.



First screen pops up very fast and disappears. It may not be visible!

Figure 6 – Options Installed Screen

12. Select [OK]

13. Select [Quit]

14. Select **[Quit]** a second time.

A window pops-up displaying: "Application should be shutdown and restarted for installed/removed options to take effect."

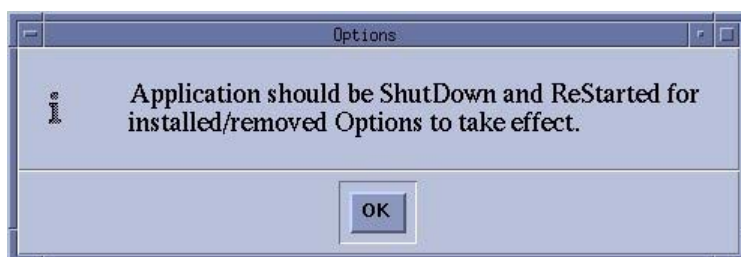


Figure 7 – Screen requesting Applications Shutdown

15. Select **[OK]**. If the install box does not go away, close the box.
16. Press the button on the DVD-RAM drive on SCSI Tower to eject the option DVD.
17. Restart applications for the options to take effect.


3. Verification of the Options Install and final checks

Verification of options install


Verify that the Neuro 3D Filter button appears in the Browser on Image Works.

- Select New Patient
- In the Protocol Number field, Type in the number 21.10 and hit Enter


See picture below

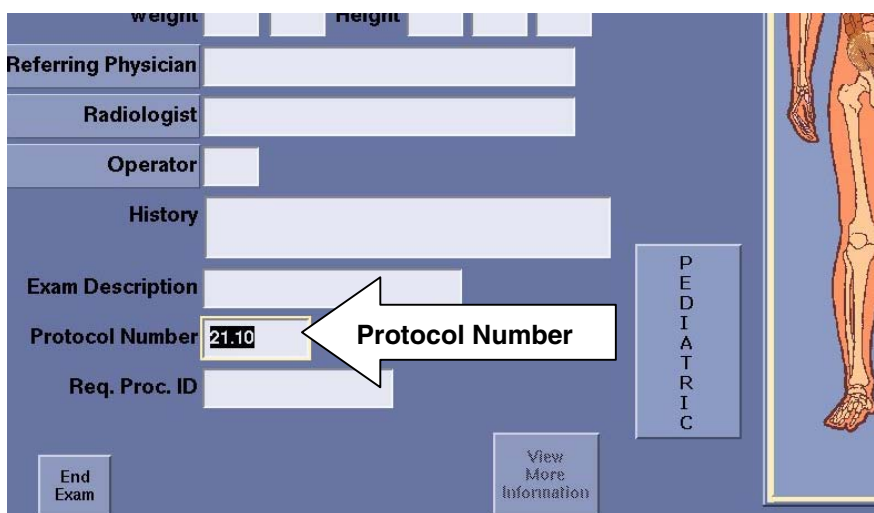


- Select Next Series

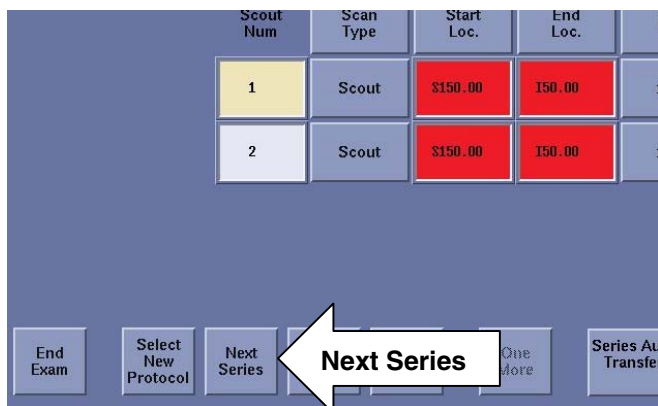

- Select Recon Tab
- Select button right below Auto Apps

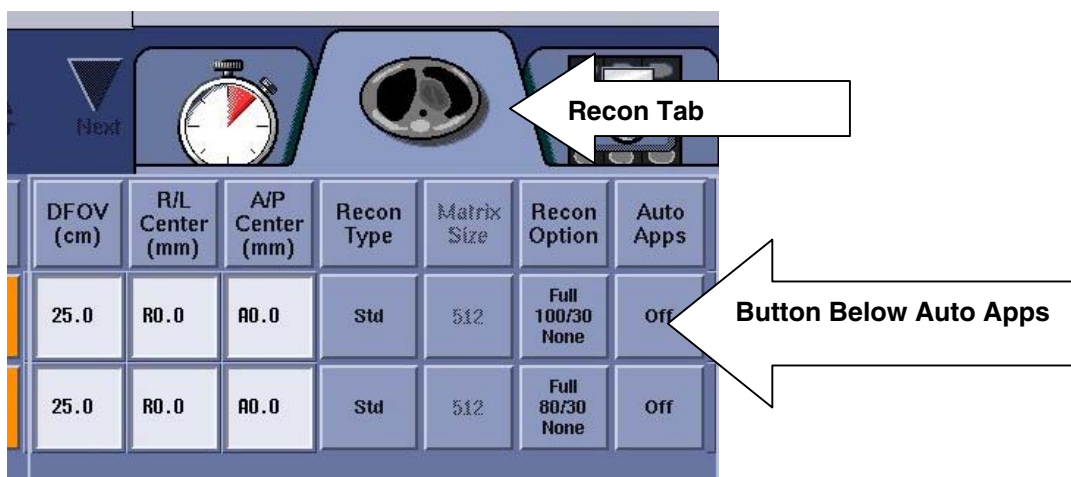

- Verify that Neuro 3D appear along with Direct 3D

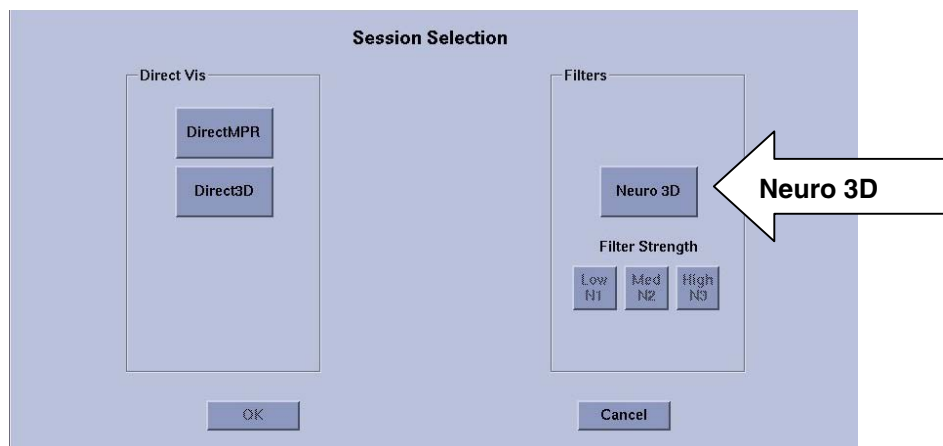

System Check:

Perform sanity scans to ensure proper system performance before returning it to the customer:

- Successfully perform a Scout scan.
- Successfully perform a Helical scan.
- Successfully perform an Axial scan

4. FMI Closure

Close FMI per regional procedure.

Remove this FMI document from the customer site.