



Technical Publication

Direction 2243318-100

Revision 08

GE Medical Systems LightSpeed 2.X - Functional Interconnects

The information in this service manual applies to the following
LightSpeed 3.X CT systems:

- LightSpeed Ultra (8-slice, MDAS)
- LightSpeed Plus (4-slice, MDAS)
- LightSpeed QX/i (4-slice, MDAS)

ADVANCED SERVICE DOCUMENTATION
PROPERTY OF GE
FOR GE SERVICE PERSONNEL ONLY
NO RIGHTS LICENSED
DO NOT USE OR COPY
DISCLOSURE TO THIRD PARTIES PROHIBITED

Copyright © 2000-2005 by General Electric Company.
All rights reserved.



PROPRIETARY TO GE MEDICAL SYSTEMS

Permission to use this Advanced Service Software and related documentation (herein called the material) by persons other than GE Medical Systems employees is provided only under an Advanced Service Package License relating specifically to this Proprietary Material. This is a different agreement from the one under which operating and basic service software is lessened. A license to use operating or basic service software does not extend to or cover this software or related documentation.

If you are a GE Medical Systems employee or a customer who has entered into such a license agreement with GE Medical Systems to use this proprietary software, you are authorized to use this Material according to the conditions stated in your license agreement. However, you do not have the permission of GE Medical Systems to alter, decompose or reverse-assemble the software, and unless you are a GE employee, you may not copy the Material. The Material is protected by Copyright and Trade Secret laws; the violation of which can result in civil damages and criminal prosecution.

If you are not party to such a license agreement or a GE Medical Systems Employee, you must exit this Material now.

IMPORTANT PRECAUTIONS

LANGUAGE

WARNING

- THIS SERVICE MANUAL IS AVAILABLE IN ENGLISH ONLY.
- IF A CUSTOMER’S SERVICE PROVIDER REQUIRES A LANGUAGE OTHER THAN ENGLISH, IT IS THE CUSTOMER’S RESPONSIBILITY TO PROVIDE TRANSLATION SERVICES.
- DO NOT ATTEMPT TO SERVICE THE EQUIPMENT UNLESS THIS SERVICE MANUAL HAS BEEN CONSULTED AND IS UNDERSTOOD.
- FAILURE TO HEED THIS WARNING MAY RESULT IN INJURY TO THE SERVICE PROVIDER, OPERATOR OR PATIENT FROM ELECTRIC SHOCK, MECHANICAL OR OTHER HAZARDS.

AVERTISSEMENT

- CE MANUEL DE MAINTENANCE N’EST DISPONIBLE QU’EN ANGLAIS.
- SI LE TECHNICIEN DU CLIENT A BESOIN DE CE MANUEL DANS UNE AUTRE LANGUE QUE L’ANGLAIS, C’EST AU CLIENT QU’IL INCOMBE DE LE FAIRE TRADUIRE.
- NE PAS TENTER D’INTERVENTION SUR LES ÉQUIPEMENTS TANT QUE LE MANUEL SERVICE N’A PAS ÉTÉ CONSULTÉ ET COMPRIS.
- LE NON-RESPECT DE CET AVERTISSEMENT PEUT ENTRAÎNER CHEZ LE TECHNICIEN, L’OPÉRATEUR OU LE PATIENT DES BLESSURES DUES À DES DANGERS ÉLECTRIQUES, MÉCANIQUES OU AUTRES.

WARNUNG

- DIESES KUNDENDIENST-HANDBUCH EXISTIERT NUR IN ENGLISCHER SPRACHE.
- FALLS EIN FREMDER KUNDENDIENST EINE ANDERE SPRACHE BENÖTIGT, IST ES AUFGABE DES KUNDEN FÜR EINE ENTSPRECHENDE ÜBERSETZUNG ZU SORGEN.
- VERSUCHEN SIE NICHT, DAS GERÄT ZU REPARIEREN, BEVOR DIESES KUNDENDIENST-HANDBUCH ZU RATE GEZOGEN UND VERSTANDEN WURDE.
- WIRD DIESE WARNUNG NICHT BEACHTET, SO KANN ES ZU VERLETZUNGEN DES KUNDENDIENSTTECHNIKERS, DES BEDIENERS ODER DES PATIENTEN DURCH ELEKTRISCHE SCHLÄGE, MECHANISCHE ODER SONSTIGE GEFAHREN KOMMEN.

AVISO

- ESTE MANUAL DE SERVICIO SÓLO EXISTE EN INGLÉS.
- SI ALGÚN PROVEEDOR DE SERVICIOS AJENO A GEMS SOLICITA UN IDIOMA QUE NO SEA EL INGLÉS, ES RESPONSABILIDAD DEL CLIENTE OFRECER UN SERVICIO DE TRADUCCIÓN.
- NO SE DEBERÁ DAR SERVICIO TÉCNICO AL EQUIPO, SIN HABER CONSULTADO Y COMPRENDIDO ESTE MANUAL DE SERVICIO.
- LA NO OBSERVANCIA DEL PRESENTE AVISO PUEDE DAR LUGAR A QUE EL PROVEEDOR DE SERVICIOS, EL OPERADOR O EL PACIENTE SUFRAN LESIONES PROVOCADAS POR CAUSAS ELÉCTRICAS, MECÁNICAS O DE OTRA NATURALEZA.

ATENÇÃO

- ESTE MANUAL DE ASSISTÊNCIA TÉCNICA SÓ SE ENCONTRA DISPONÍVEL EM INGLÊS.
- SE QUALQUER OUTRO SERVIÇO DE ASSISTÊNCIA TÉCNICA, QUE NÃO A GEMS, SOLICITAR ESTES MANUAIS NOUTRO IDIOMA, É DA RESPONSABILIDADE DO CLIENTE FORNECER OS SERVIÇOS DE TRADUÇÃO.
- NÃO TENDE REPARAR O EQUIPAMENTO SEM TER CONSULTADO E COMPREENDIDO ESTE MANUAL DE ASSISTÊNCIA TÉCNICA.
- O NÃO CUMPRIMENTO DESTA AVISO PODE POR EM PERIGO A SEGURANÇA DO TÉCNICO, OPERADOR OU PACIENTE DEVIDO A‘ CHOQUES ELÉTRICOS, MECÂNICOS OU OUTROS.

AVVERTENZA

- IL PRESENTE MANUALE DI MANUTENZIONE È DISPONIBILE SOLTANTO IN INGLESE.
- SE UN ADDETTO ALLA MANUTENZIONE ESTERNO ALLA GEMS RICHIEDE IL MANUALE IN UNA LINGUA DIVERSA, IL CLIENTE È TENUTO A PROVVEDERE DIRETTAMENTE ALLA TRADUZIONE.
- SI PROCEDA ALLA MANUTENZIONE DELL’APPARECCHIATURA SOLO DOPO AVER CONSULTATO IL PRESENTE MANUALE ED AVERNE COMPRESO IL CONTENUTO.
- NON TENERE CONTO DELLA PRESENTE AVVERTENZA POTREBBE FAR COMPIERE OPERAZIONI DA CUI DERIVINO LESIONI ALL’ADDETTO ALLA MANUTENZIONE, ALL’UTILIZZATORE ED AL PAZIENTE PER FOLGORAZIONE ELETTRICA, PER URTI MECCANICI OD ALTRI RISCHI.

警告

このサービスマニュアルには英語版しかありません。
GEMS以外でサービスを担当される業者が英語以外の言語を要求される場合、翻訳作業はその業者の責任で行うものとさせていただきます。
このサービスマニュアルを熟読し理解せずに、装置のサービスを行わないで下さい。
この警告に従わない場合、サービスを担当される方、操作員あるいは患者さんが、感電や機械的又はその他の危険により負傷する可能性があります。

注意:

本维修手册仅存有英文本。
非 GEMS 公司的维修员要求非英文本的维修手册时，客户需自行负责翻译。
未详细阅读和完全了解本手册之前，不得进行维修。
忽略本注意事项会对维修员，操作员或病人造成触电，机械伤害或其他伤害。

DAMAGE IN TRANSPORTATION

All packages should be closely examined at time of delivery. If damage is apparent, have notation “damage in shipment” written on all copies of the freight or express bill before delivery is accepted or “signed for” by a General Electric representative or a hospital receiving agent. Whether noted or concealed, damage MUST be reported to the carrier immediately upon discovery, or in any event, within 14 days after receipt, and the contents and containers held for inspection by the carrier. A transportation company will not pay a claim for damage if an inspection is not requested within this 14 day period.

To file a report:

- Call 1-800-548-3366 and use option 8.
- Fill out a report on <http://us44hdd21/sctq/InstallFulfill/InstalFulfillment.htm>
- Contact your local service coordinator for more information on this process.

Rev. Jan. 5, 2005

CERTIFIED ELECTRICAL CONTRACTOR STATEMENT

All electrical Installations that are preliminary to positioning of the equipment at the site prepared for the equipment shall be performed by licensed electrical contractors. In addition, electrical feeds into the Power Distribution Unit shall be performed by licensed electrical contractors. Other connections between pieces of electrical equipment, calibrations and testing shall be performed by qualified GE Medical personnel. The products involved (and the accompanying electrical installations) are highly sophisticated, and special engineering competence is required. In performing all electrical work on these products, GE will use its own specially trained field engineers. All of GE’s electrical work on these products will comply with the requirements of the applicable electrical codes.

The purchaser of GE equipment shall only utilize qualified personnel (i.e., GE’s field engineers, personnel of third-party service companies with equivalent training, or licensed electricians) to perform electrical servicing on the equipment.

IMPORTANT...X-RAY PROTECTION

X-ray equipment if not properly used may cause injury. Accordingly, the instructions herein contained should be thoroughly read and understood by everyone who will use the equipment before you attempt to place this equipment in operation. The General Electric Company, Medical Systems Group, will be glad to assist and cooperate in placing this equipment in use.

Although this apparatus incorporates a high degree of protection against x-radiation other than the useful beam, no practical design of equipment can provide complete protection. Nor can any practical design compel the operator to take adequate precautions to prevent the possibility of any persons carelessly exposing themselves or others to radiation.

It is important that anyone having anything to do with x-radiation be properly trained and fully acquainted with the recommendations of the National Council on Radiation Protection and Measurements as published in NCRP Reports available from NCRP Publications, 7910 Woodmont Avenue, Room 1016, Bethesda, Maryland 20814, and of the International Commission on Radiation Protection, and take adequate steps to protect against injury.

The equipment is sold with the understanding that the General Electric Company, Medical Systems Group, its agents, and representatives have no responsibility for injury or damage which may result from improper use of the equipment.

Various protective materials and devices are available. It is urged that such materials or devices be used.

LITHIUM BATTERY CAUTIONARY STATEMENTS



CAUTION
Risk of
Explosion

Danger of explosion if battery is incorrectly replaced. Replace only with the same or equivalent type recommended by the manufacturer. Discard used batteries according to the manufacturer’s instructions.



ATTENTION
Danger
d’Explosion

Il y a danger d’explosion s’il y a remplacement incorrect de la batterie. Remplacer uniquement avec une batterie du même type ou d’un type recommandé par le constructeur. Mettre au rebut les batteries usagées conformément aux instructions du fabricant.

OMISSIONS & ERRORS

Customers, please contact your GE Sales or Service representatives.

GE personnel, please use the GEMS PQR Process to report all omissions, errors, and defects in this publication.

Revision History

Rev	Date	Reason for change
0	11/20/00	Initial Release for LightSpeed Plus
1	02/28/01	CQA 1011476: Updated 2267773SCH to Rev 1 (Removed Elev & Tilt Drive Contactor, per PCN 220017) Added new Image Generation Functional Interconnects.
2	04/20/01	Added Interconnects for Gantry Options (2286501) & Smartview Controls (2286504)
3	07/20/01	Fixed various drawings that were missing information
4	02/12/2002	CQA 1020694: Added drawing 2151953 Rev. 0 Sheets 1 and 2.
5	06/28/2002	CQA 1020172: Added Chapter 3, including 2302517SCH
6	12/18/2003	CQA 1033884: Added drawing 2126249 to Chapter 4
7	08/17/2004	PQR 13012917: Added drawing 2123540 to Chapter 4
8	05/10/2005	PQR 13024644: Updated drawing 2267770SCH, Added drawing 2267713SCH Reformatted manual for 11x17 printing.

Table of Contents

Chapter 1	
System Communications.....	7
Section 1.0 — Gantry Slip Ring Communications - 2267732	8
Section 2.0 — Rotating Controller Interface Bus (RCIB) - 2267733.....	9
Section 3.0 — Final Scan Control / Intercom / AutoVoice - 2151953.....	10
 Chapter 2	
X-Ray Generation	13
Section 1.0 — Rotor Control - 2267317	14
Section 2.0 — Exposure Enable - 2267735	15
Section 3.0 — Rotating HVON Sense - 2267736	16
Section 4.0 — Filament Select and mA Loop - 2267737	17
Section 5.0 — Tube ID, Cooling and Pressure Sense - 2267738.....	18
Section 6.0 — Back-Up Timer - 2267739	19
Section 7.0 — kV Loop - 2267740	20
Section 8.0 — BU Contactor / Exposure Interlocks - 2267741	22
Section 9.0 — X-Ray Collimation and Filtration - 2267744	24
Section 10.0 — Final Exposure Command - 2267746	25
 Chapter 3	
Data Acquisition	27
Section 1.0 — Scan Data Trigger Generation - 2126245.....	28
Section 2.0 — Scan Data Acquisition - 2269964.....	29
Section 3.0 — Detector Heater - 2236260.....	30
Section 4.0 — DAS - 2302517	31
 Chapter 4	
Patient Positioning	33
Section 1.0 — Patient Alignment Lights - 2126249	34
Section 2.0 — Table Gantry Sync - 2267745	35
Section 3.0 — Table Elevation Drive - 2267754	36
Section 4.0 — Gantry Tilt Function - 2267756	37
Section 5.0 — Table Cradle Drive (Vertical & Longitudinal Drives) - 2123540 .	38

Chapter 5	
Axial Control	41
Section 1.0 — Axial Servo Control Loop - 2267757.....	42
Section 2.0 — Axial Drive Power Control - 2267784.....	43
 Chapter 6	
Operator I/O.....	45
Section 1.0 — Gantry Operator Interface - 2267760	46
Section 2.0 — Site X-Ray ON Light - 2267761	47
Section 3.0 — Gantry Options - 2286501.....	48
Section 4.0 — Smartview Controls - 2286504	49
Section 5.0 — RX Tilt & Gantry Console Reset - 2267765	50
 Chapter 7	
System Monitoring	51
Section 1.0 — Mains Undervoltage - 2267731.....	52
Section 2.0 — Gantry Temp. Sensors & Power Supply Monitors - 2267766....	53
 Chapter 8	
System Power Control	55
Section 1.0 — HVDC Bus Control - 2267769	56
Section 2.0 — Gantry 120VAC Control - 2267770.....	58
Section 3.0 — E-Stop Drives Reset - 2267773	59
Section 4.0 — Gantry Stationary Power and Ground - 2267713	60
 Chapter 9	
Image Generation	65
Section 1.0 — Scan Data Save	66
Section 2.0 — Data Restore.....	67
Section 3.0 — Pre-Processing, Filtered Back Projection & Post Processing..	68
Section 4.0 — Image Transfer	69

NOTES:

Chapter 1

System Communications

Section 1.0 — Gantry Slip Ring Communications - 2267732

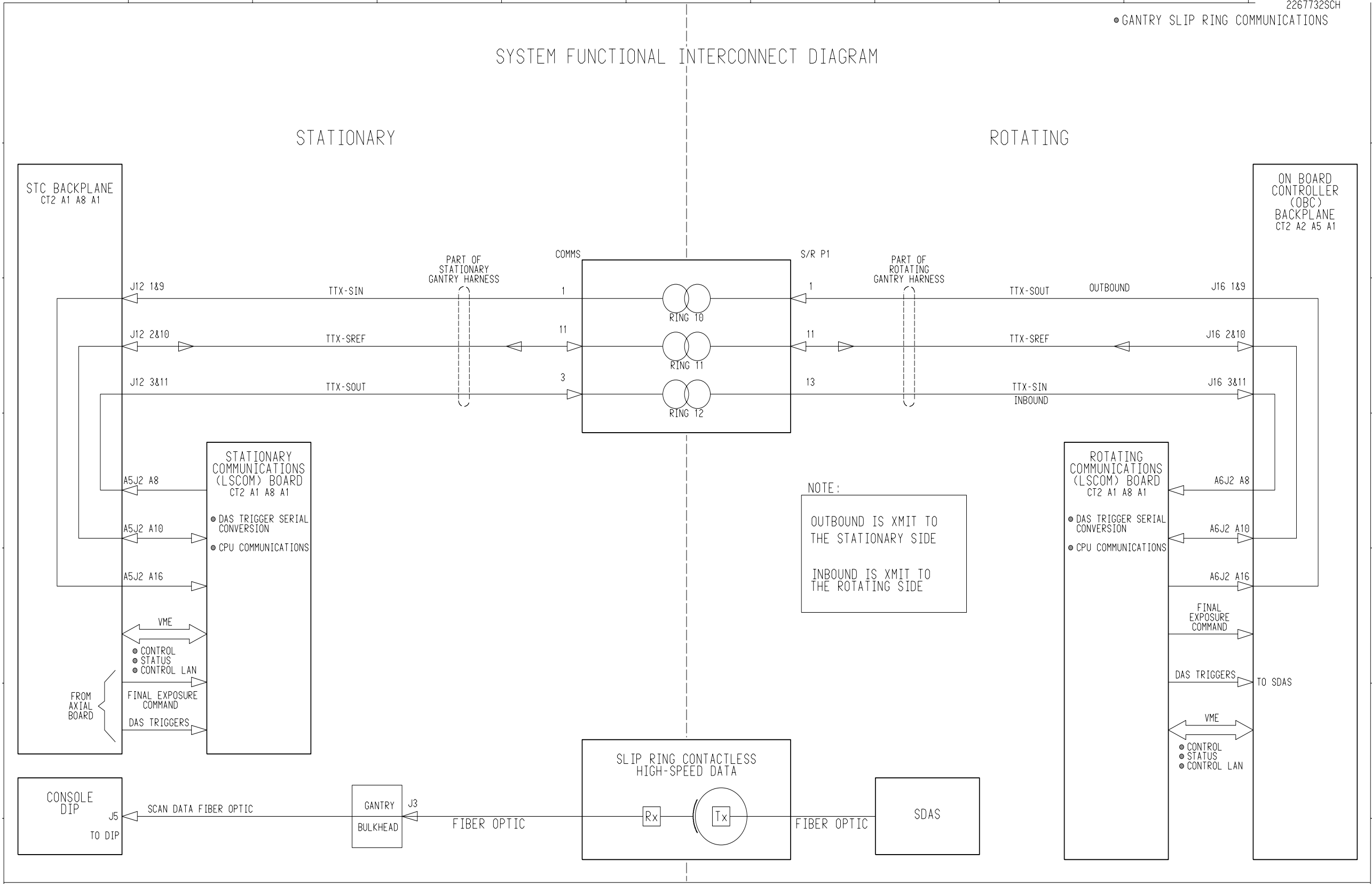


Figure 1-1 2267732 Sheet 1, Rev. 0

Section 3.0 — Final Scan Control / Intercom / AutoVoice - 2151953

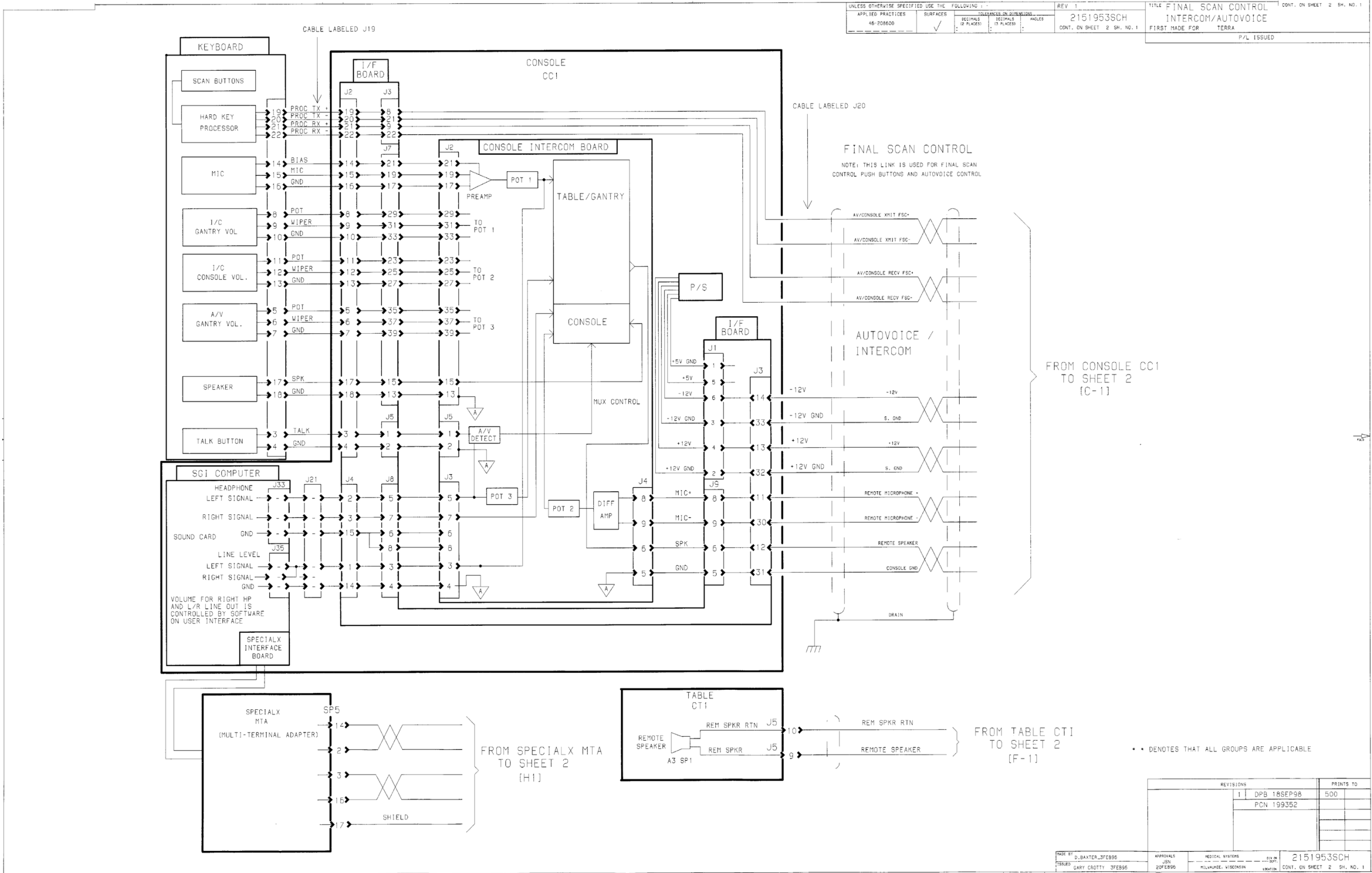


Figure 1-3 2151953 Sheet 1, Rev 1

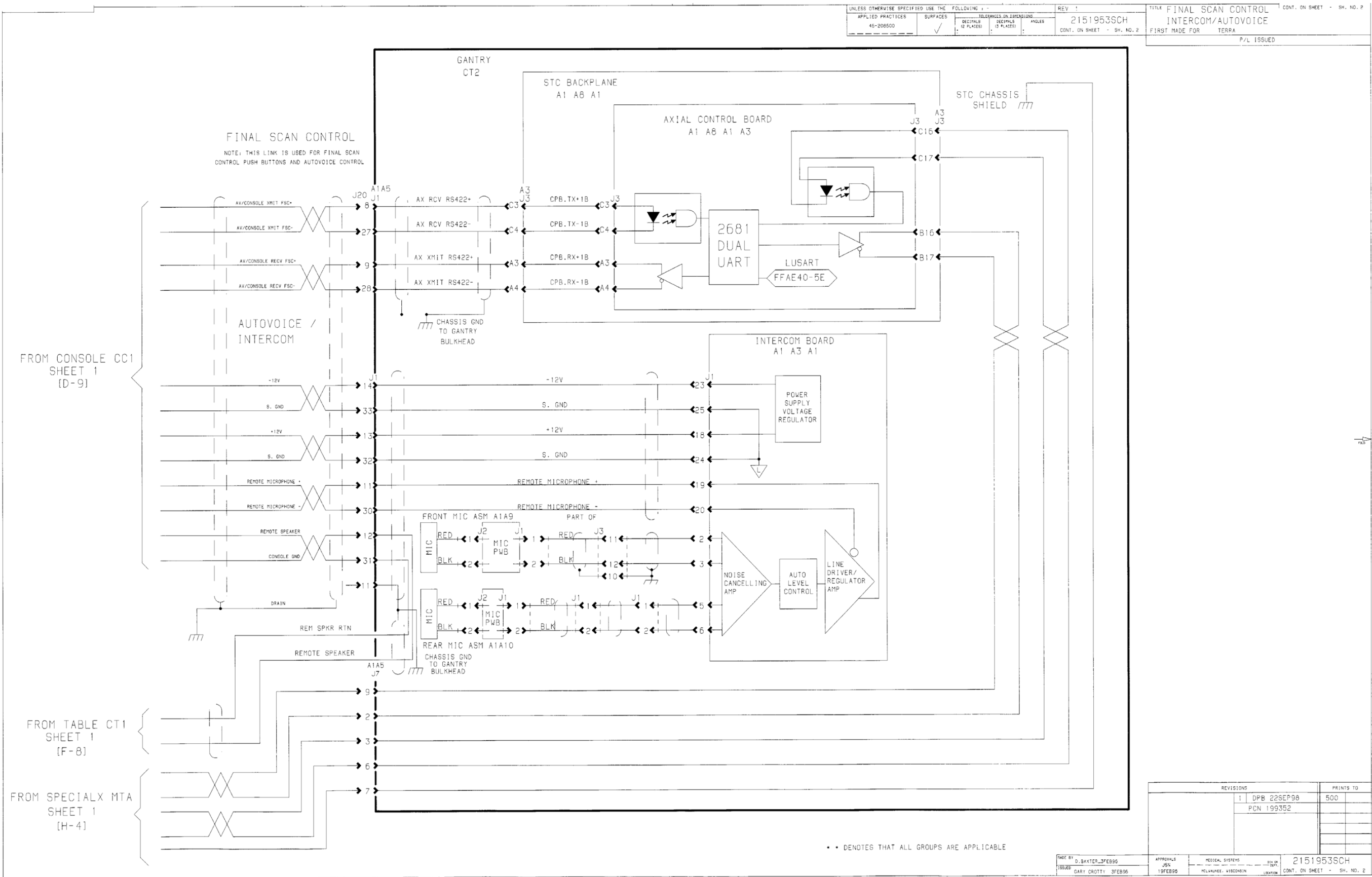


Figure 1-4 2151953 Sheet 2, Rev 1

NOTES:

Chapter 2

X-Ray Generation

Section 1.0 — Rotor Control - 2267317

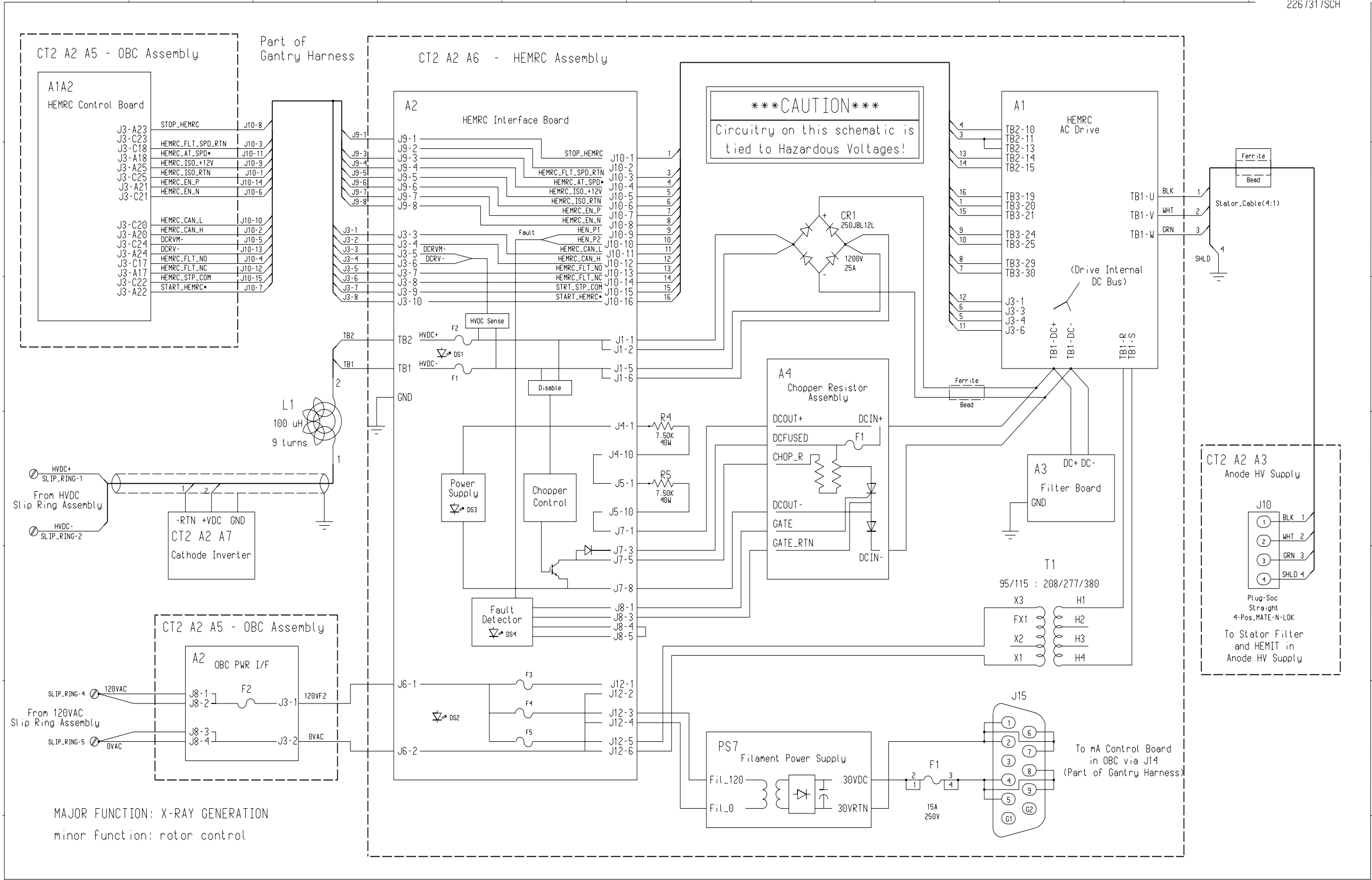


Figure 2-1 2267317 Sheet 1, Rev. 0

Section 2.0 — Exposure Enable - 2267735

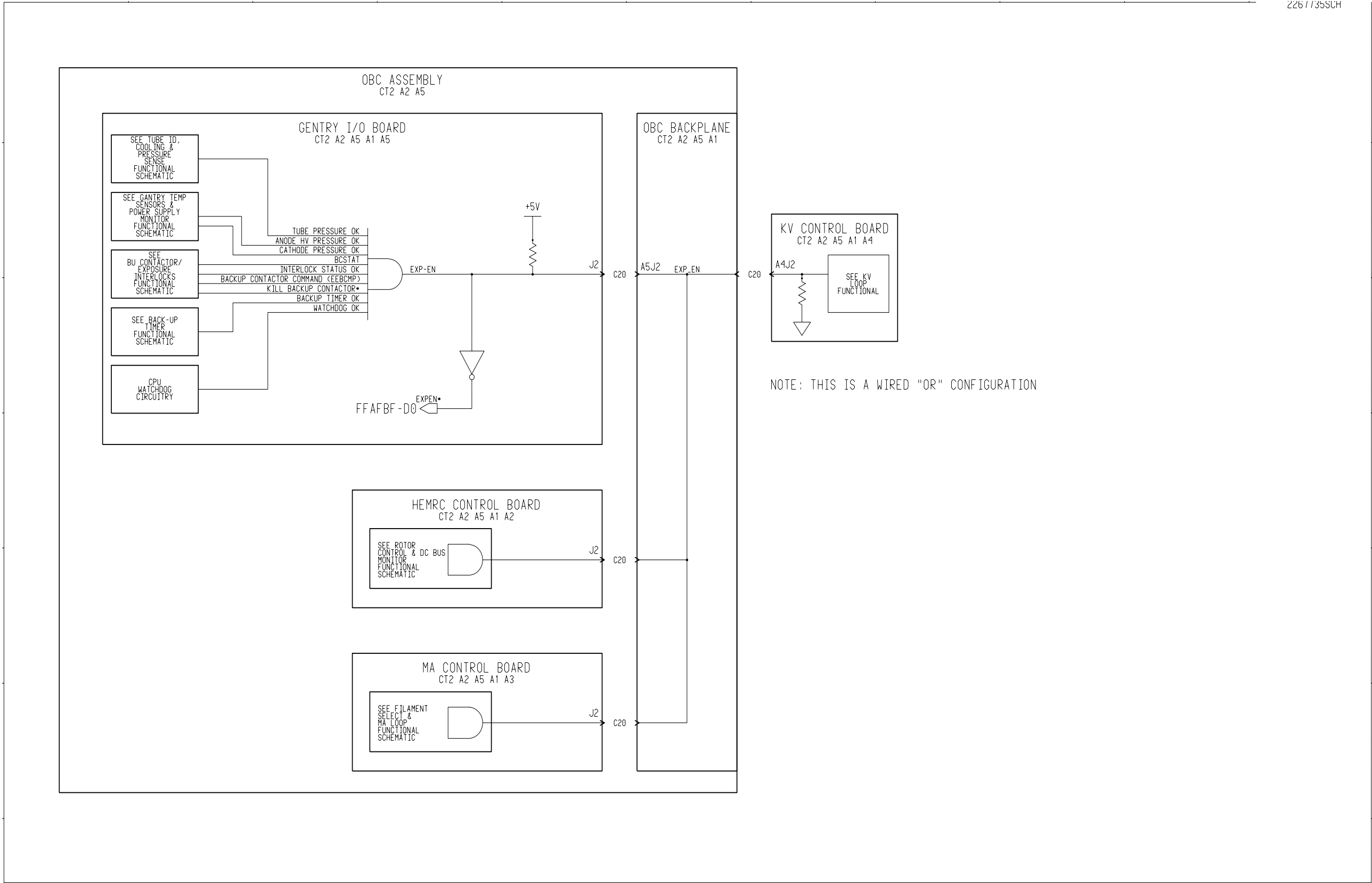


Figure 2-2 2267735 Sheet 1, Rev. 0

Section 3.0 — Rotating HVON Sense - 2267736

2267736SCH

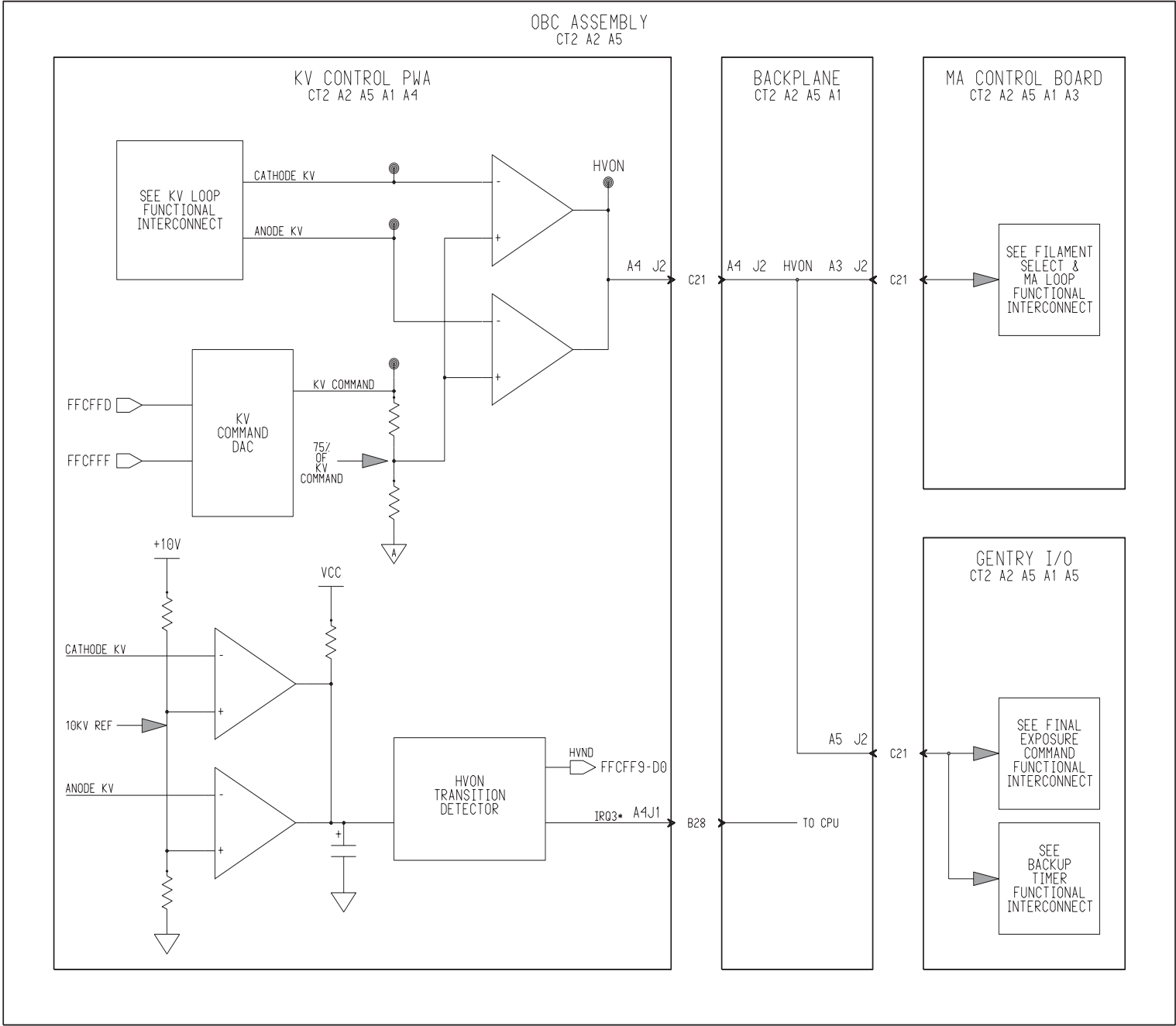


Figure 2-3 2267736 Sheet 1, Rev. 0

Section 4.0 — Filament Select and mA Loop - 2267737

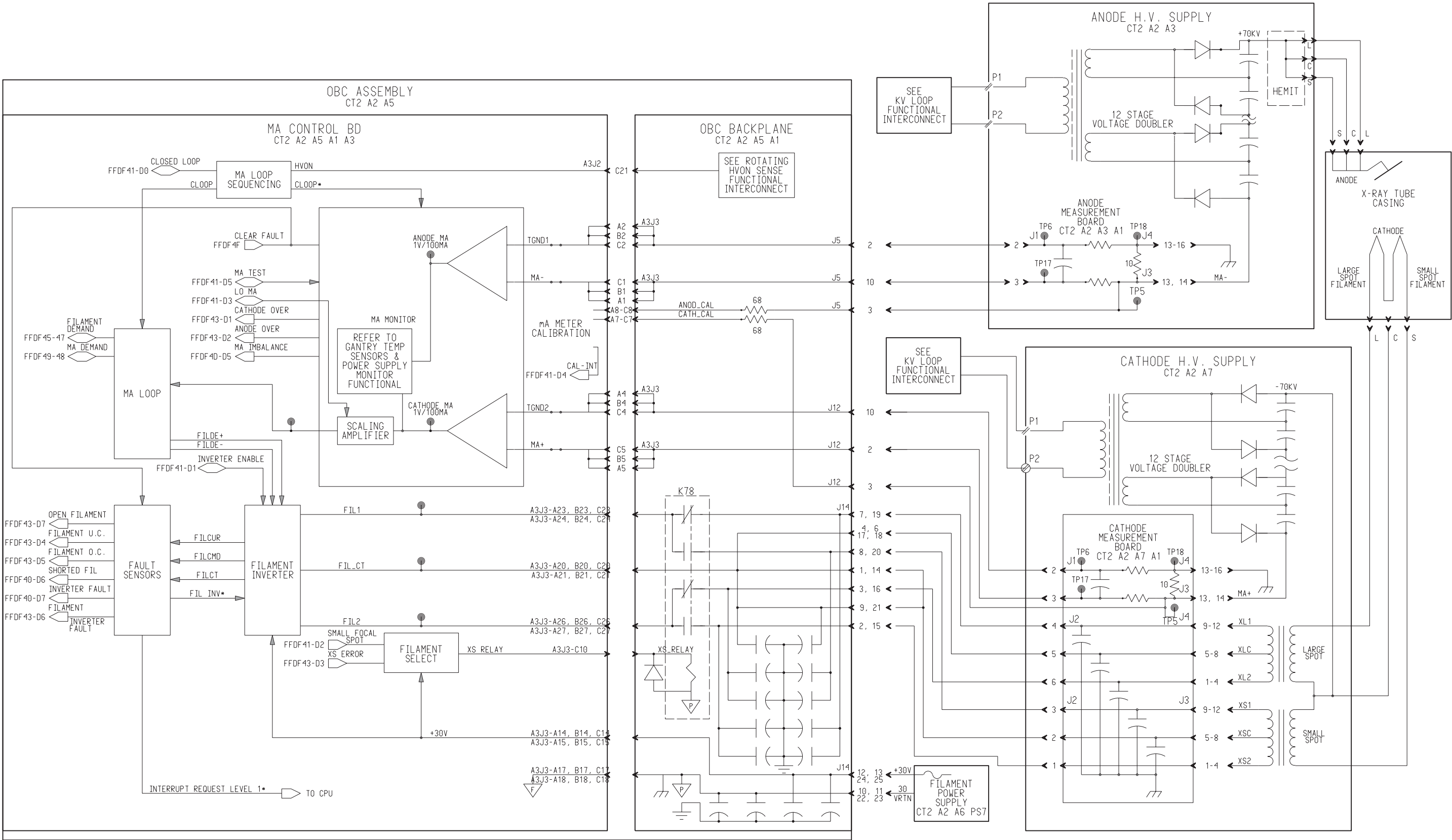


Figure 2-4 2267737 Sheet 1, Rev. 0

Section 5.0 — Tube ID, Cooling and Pressure Sense - 2267738

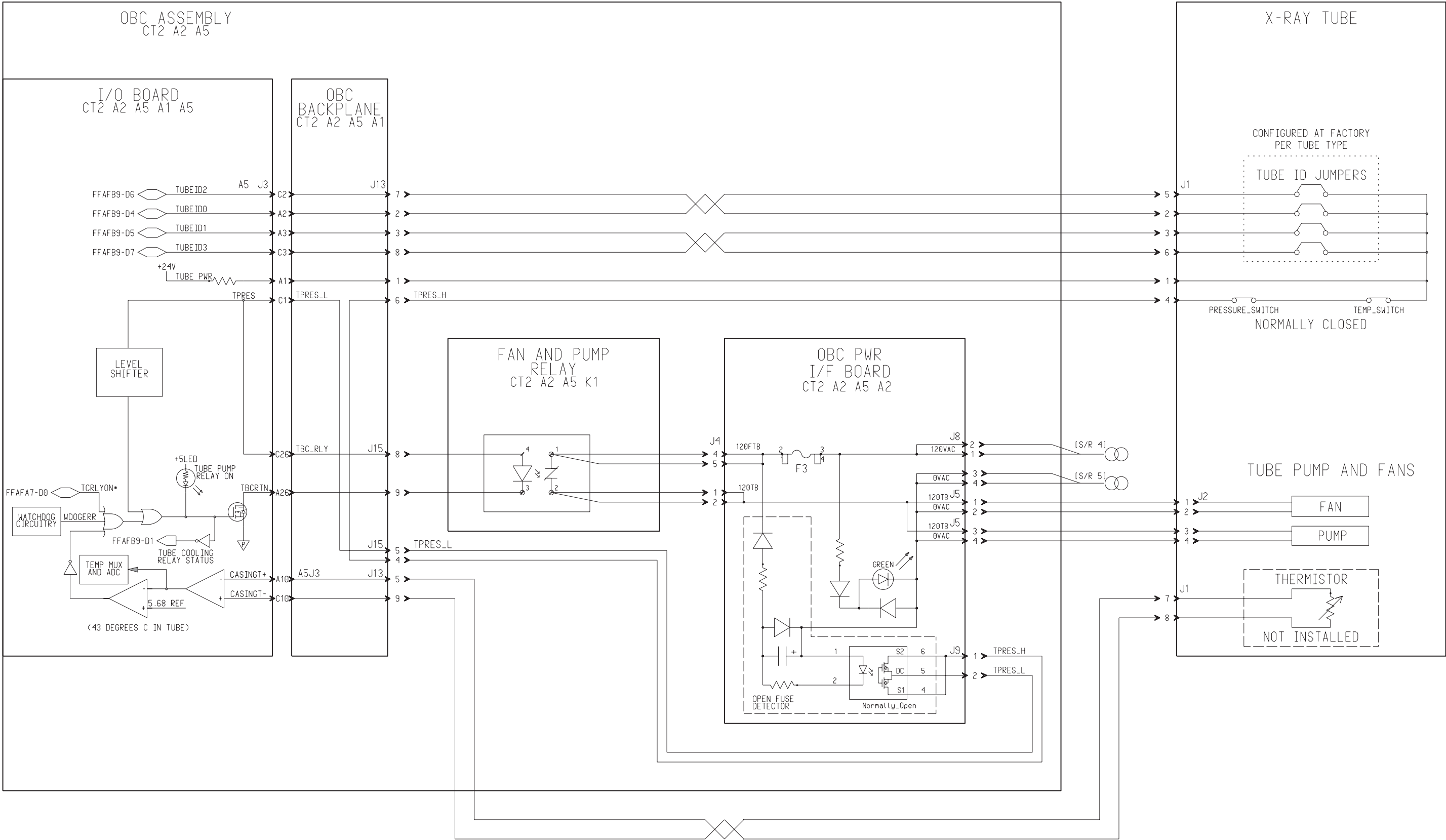


Figure 2-5 2267738 Sheet 1, Rev. 0

Section 6.0 — Back-Up Timer - 2267739


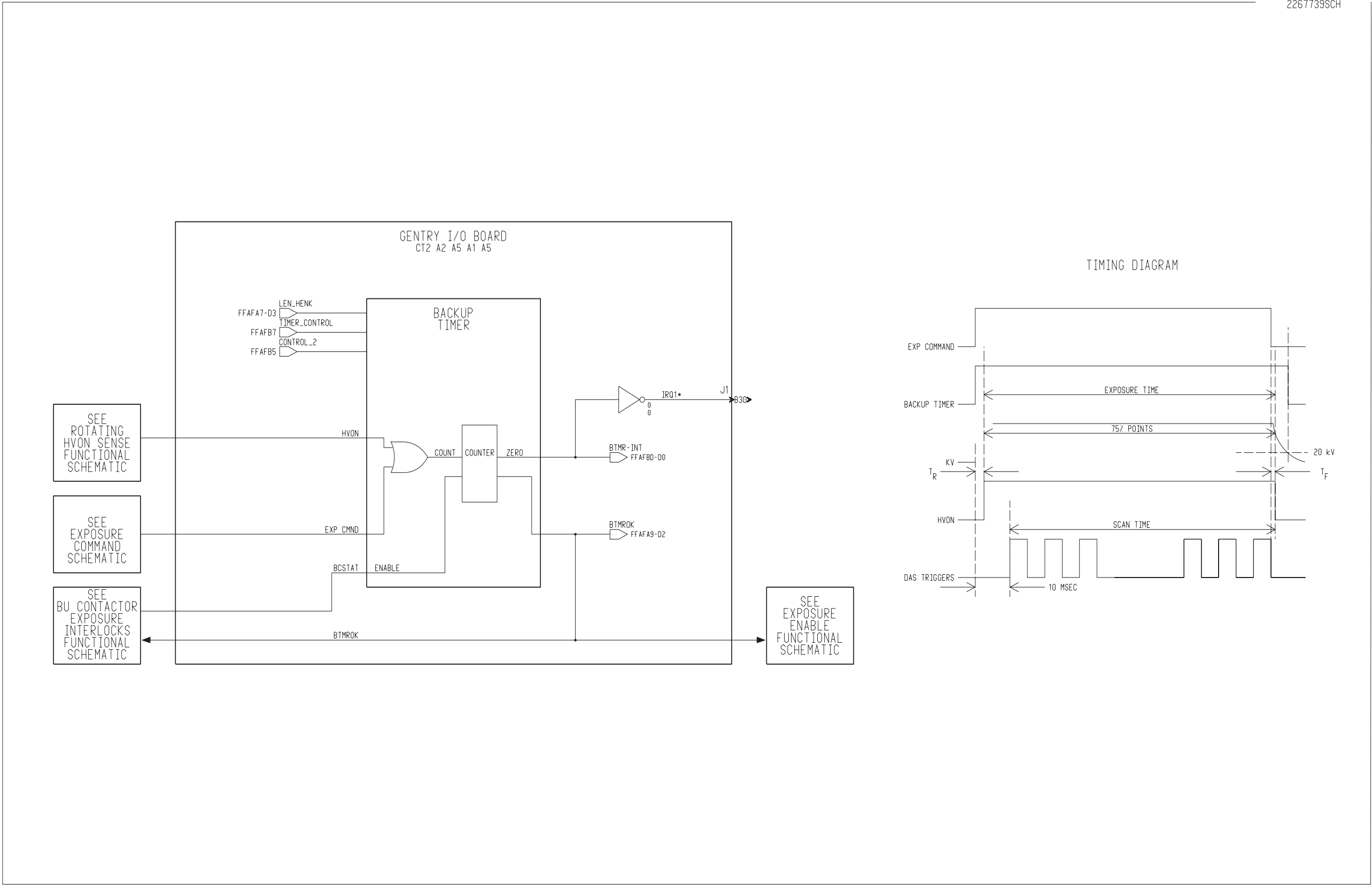


Figure 2-6 2267739 Sheet 1, Rev. 0

Section 7.0 — kV Loop - 2267740



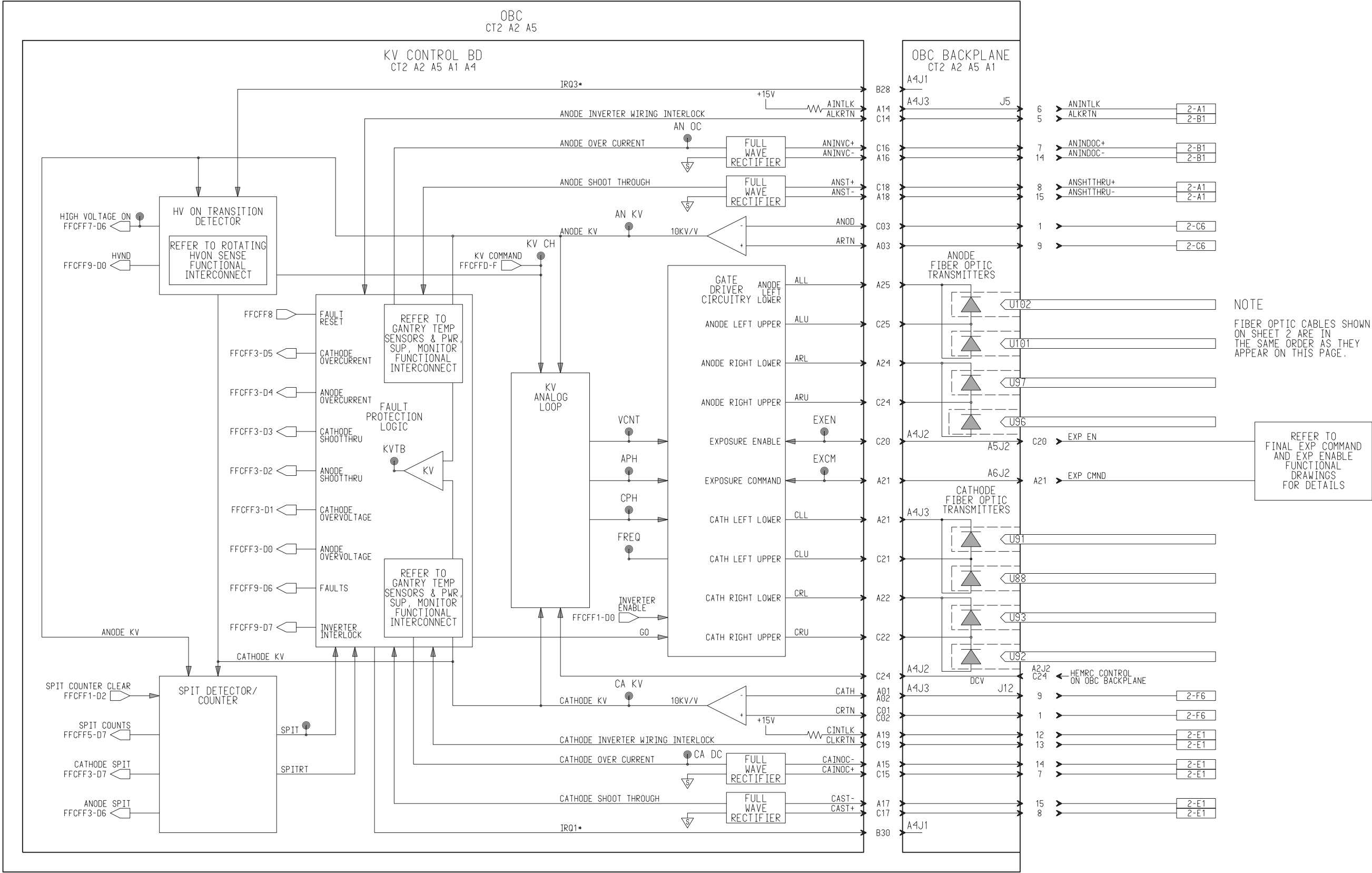

Section 8.0 — BU Contactor / Exposure Interlocks - 2267741

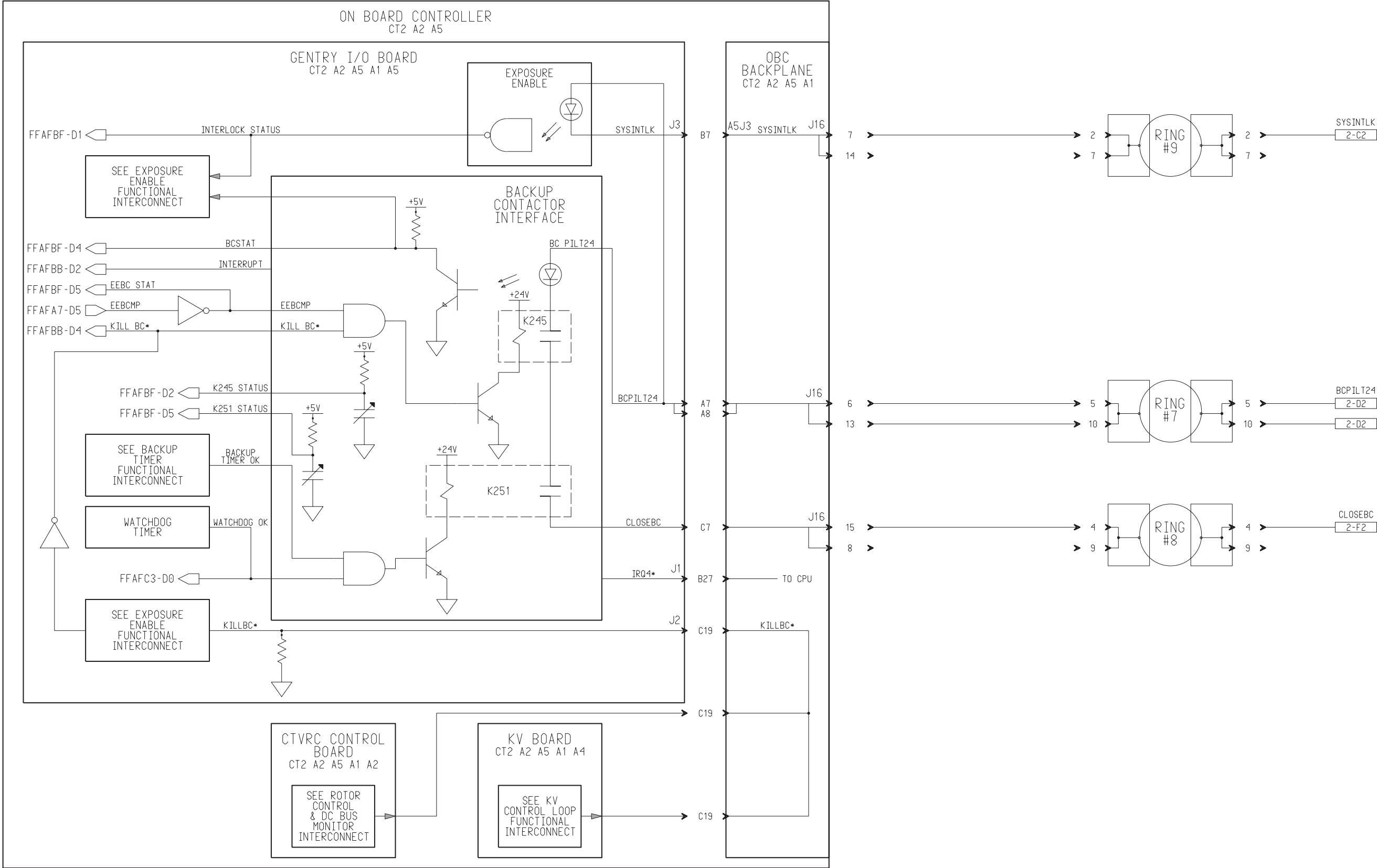


Figure 2-9 2267741 Sheet 1, Rev. 0

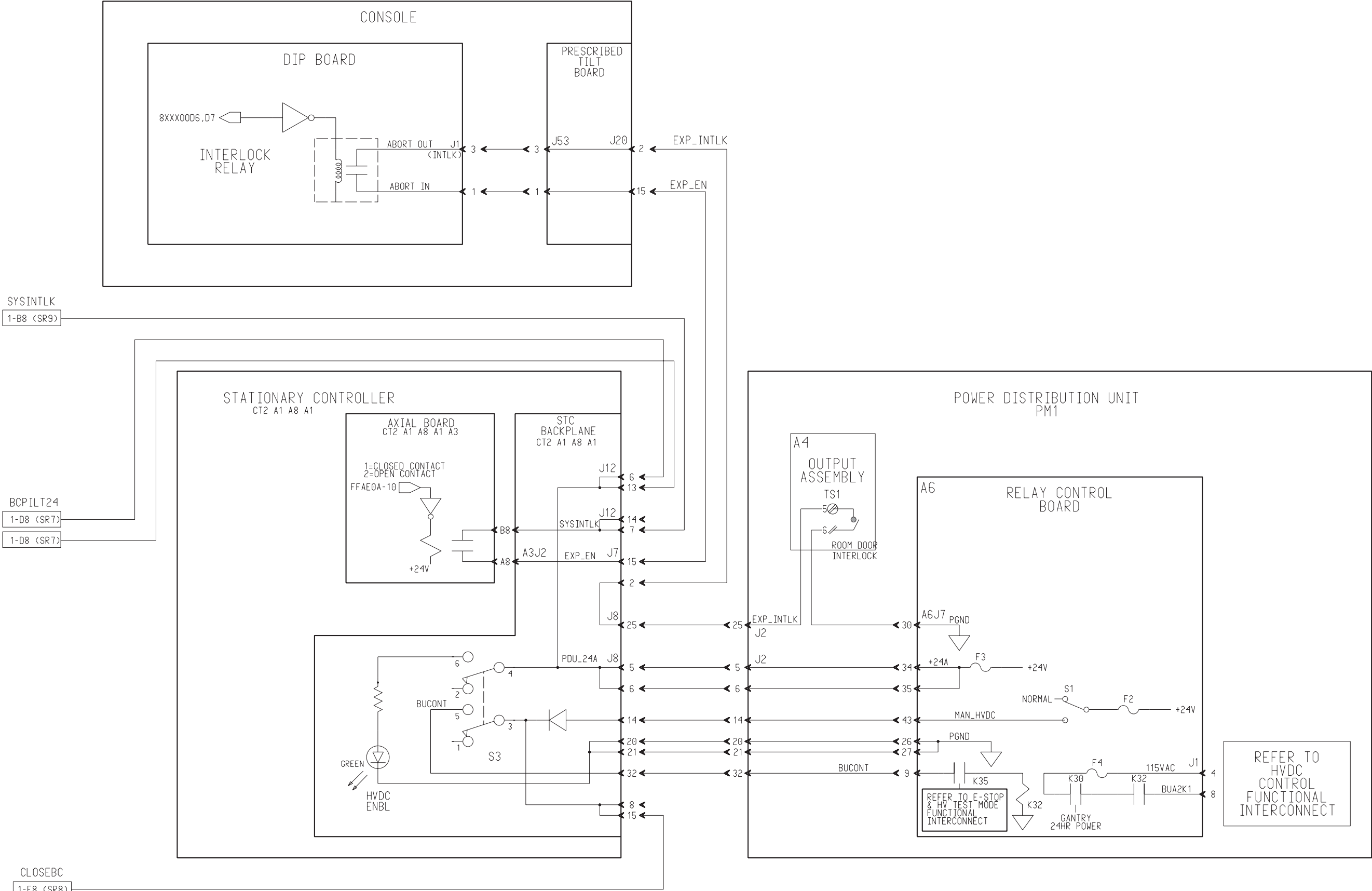


Figure 2-10 2267741 Sheet 2, Rev. 0

Section 9.0 — X-Ray Collimation and Filtration - 2267744

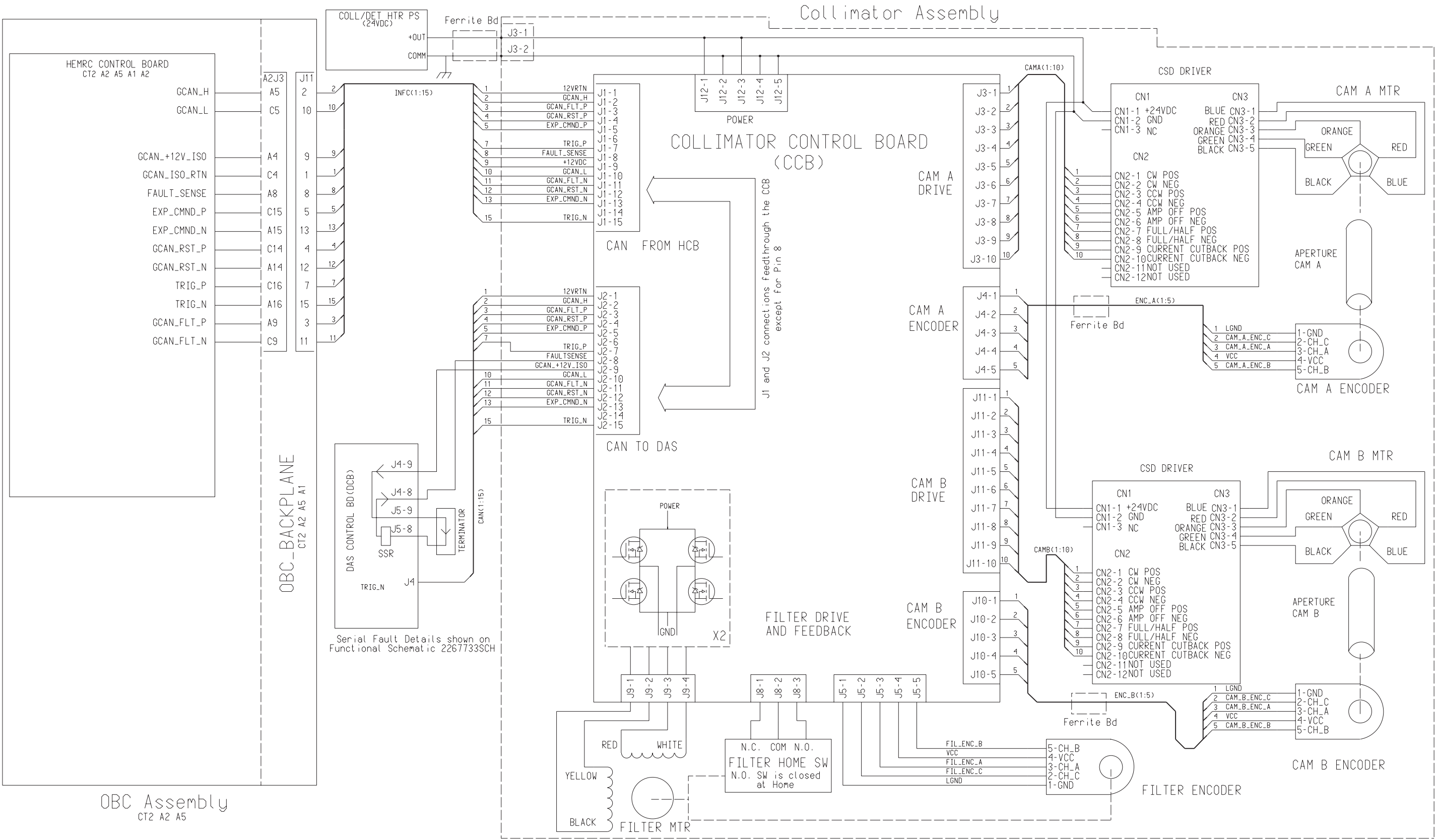


Figure 2-11 2267744 Sheet 1, Rev. 0

Section 10.0 — Final Exposure Command - 2267746

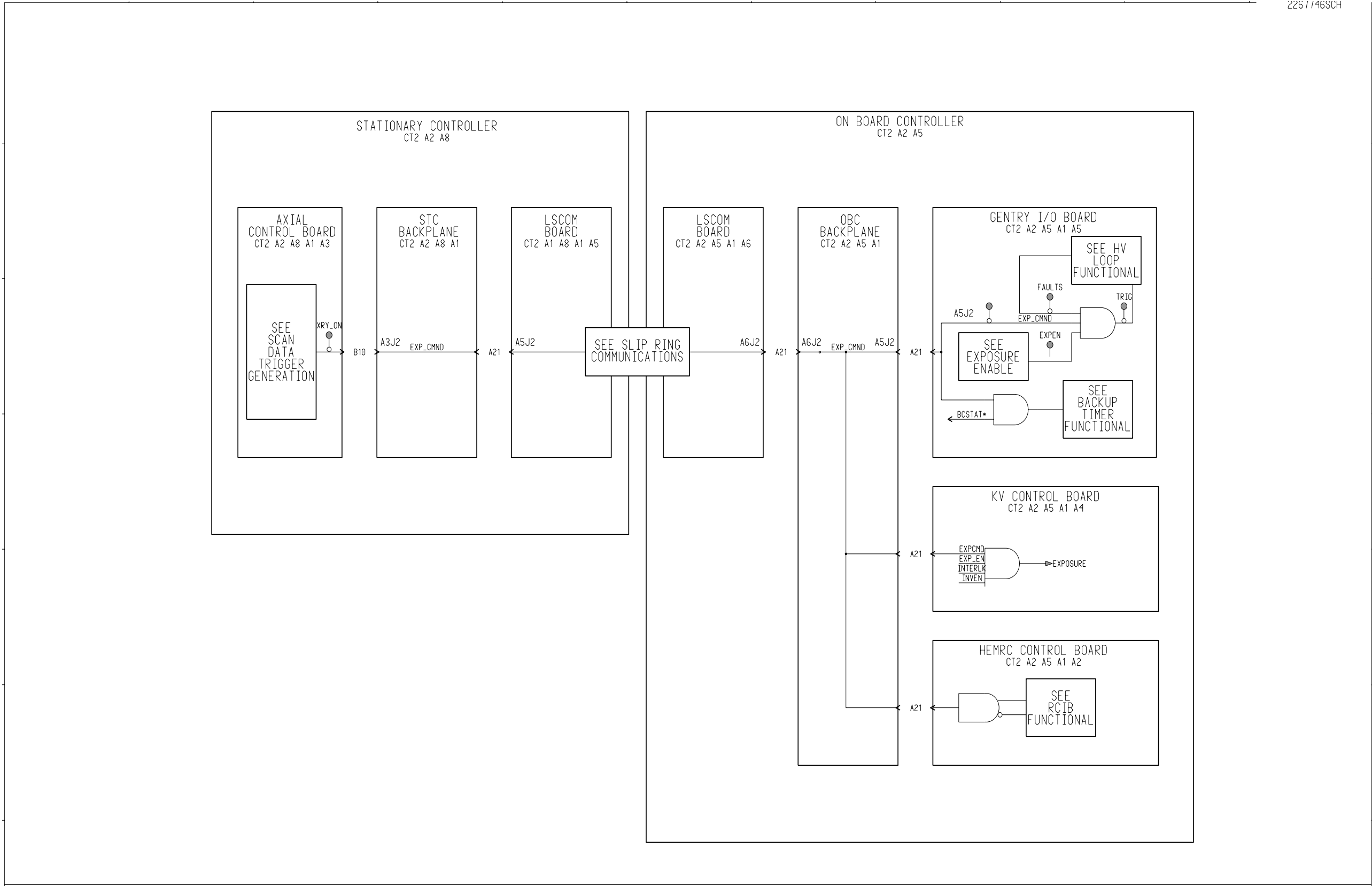


Figure 2-12 2267746 Sheet 1, Rev. 0

NOTES:

Chapter 3

Data Acquisition

Section 1.0 — Scan Data Trigger Generation - 2126245

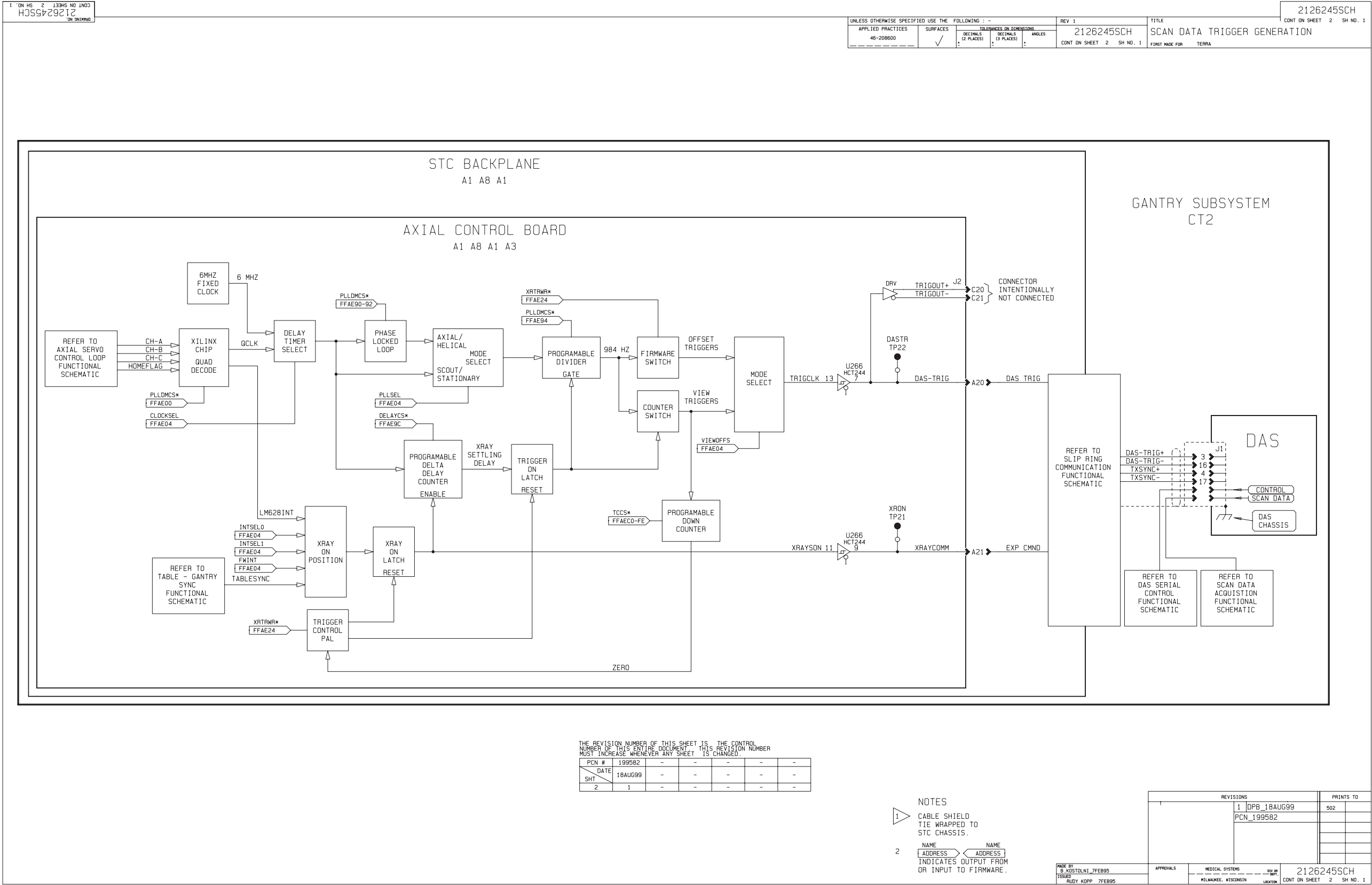


Figure 3-1 2126245 Sheet 1, Rev. 1

Section 2.0 — Scan Data Acquisition - 2269964

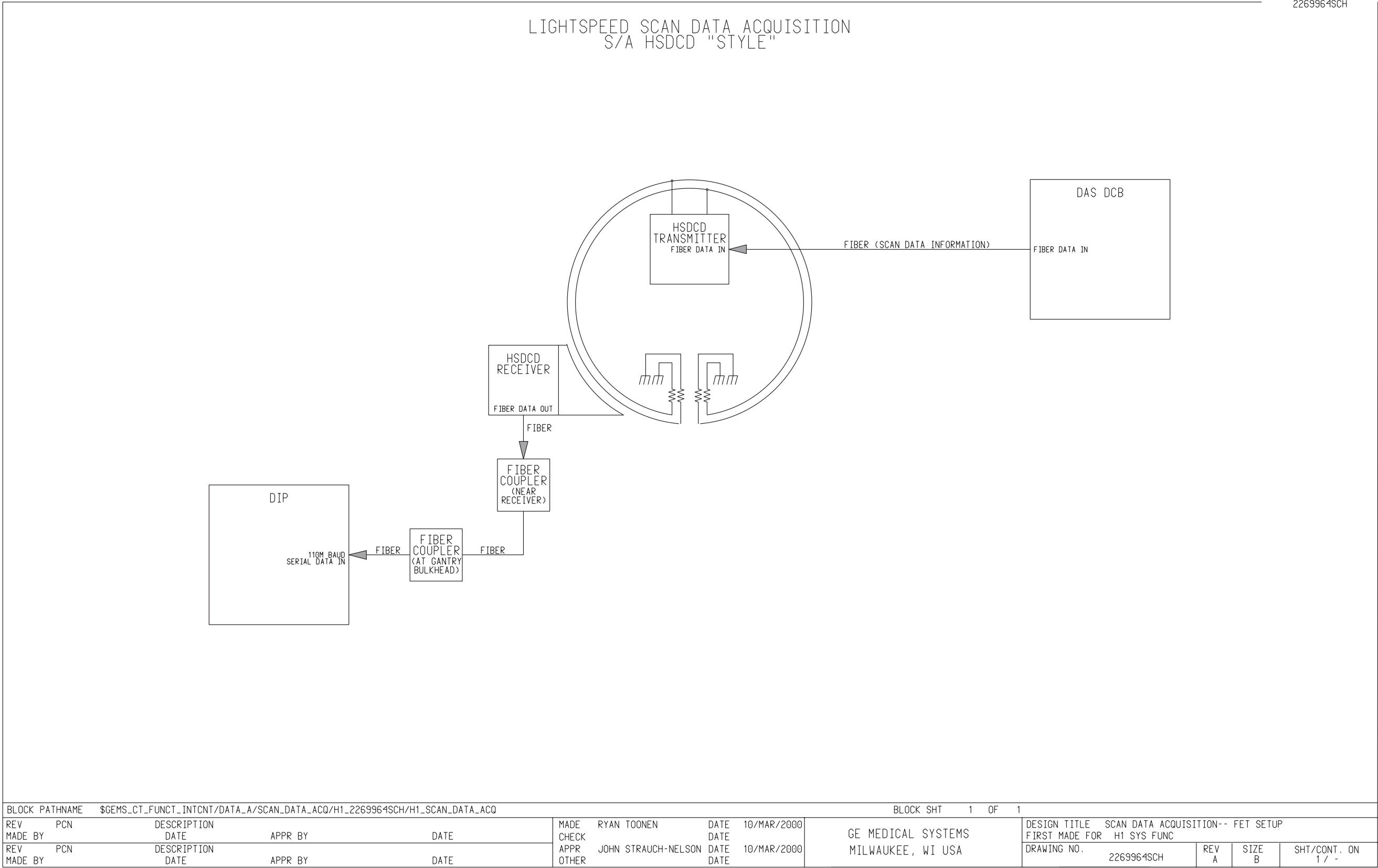


Figure 3-2 2269964 Sheet 1, Rev. A

Section 4.0 — DAS - 2302517

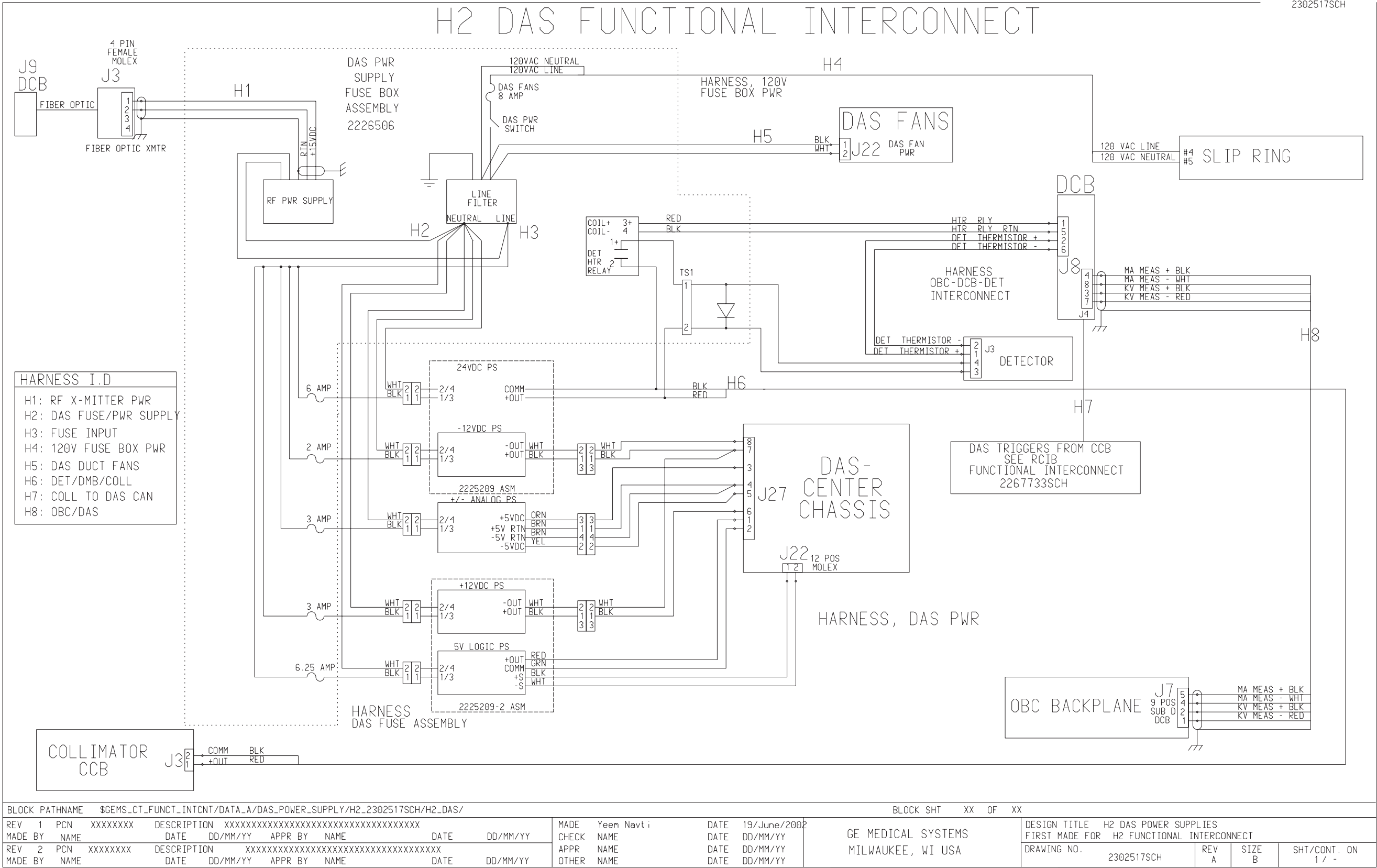



Figure 3-4 2302517 Sheet 1, Rev. A

NOTES:

Chapter 4

Patient Positioning

Section 1.0 — Patient Alignment Lights - 2126249

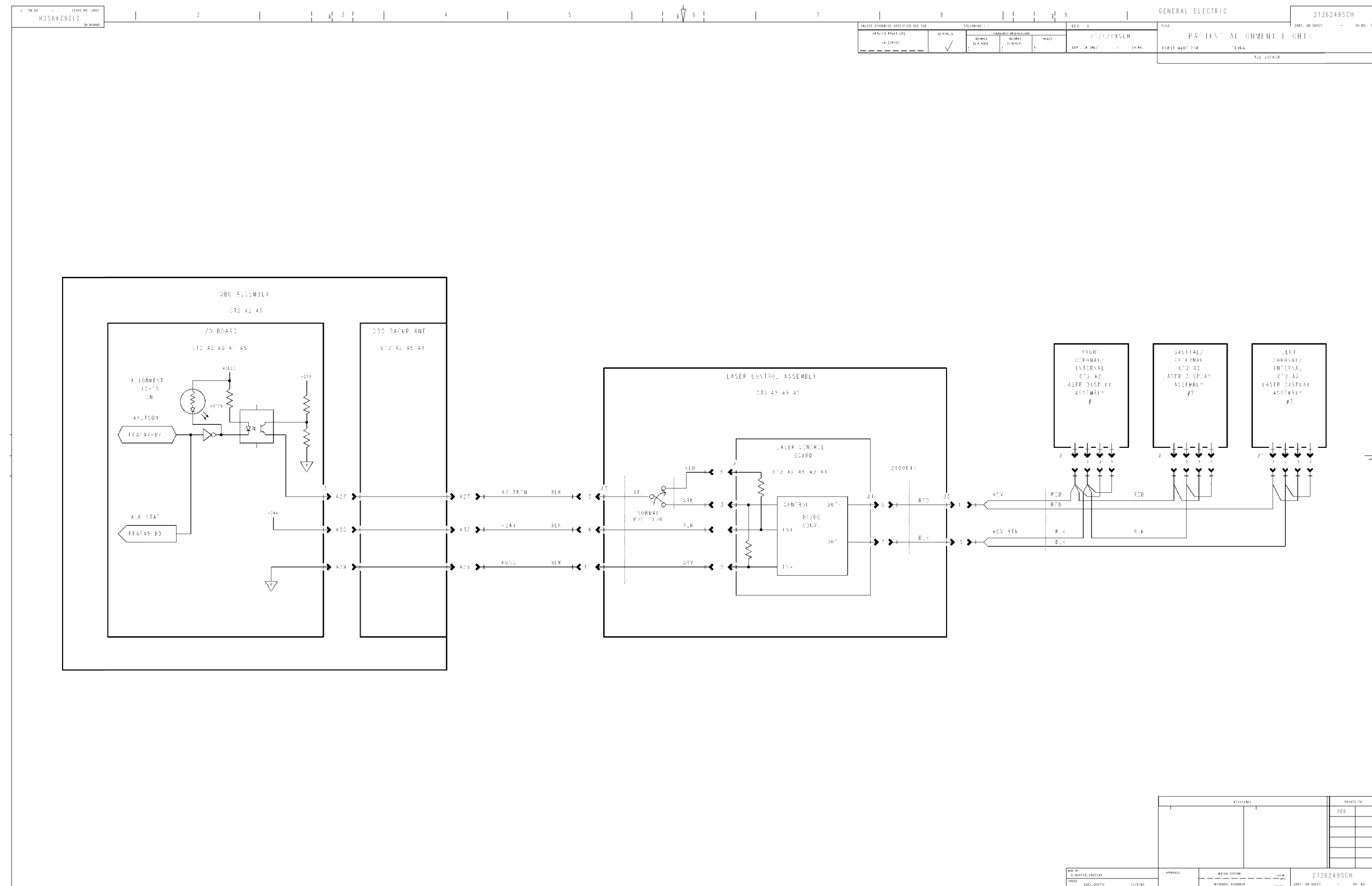

Section 2.0 — Table Gantry Sync - 2267745



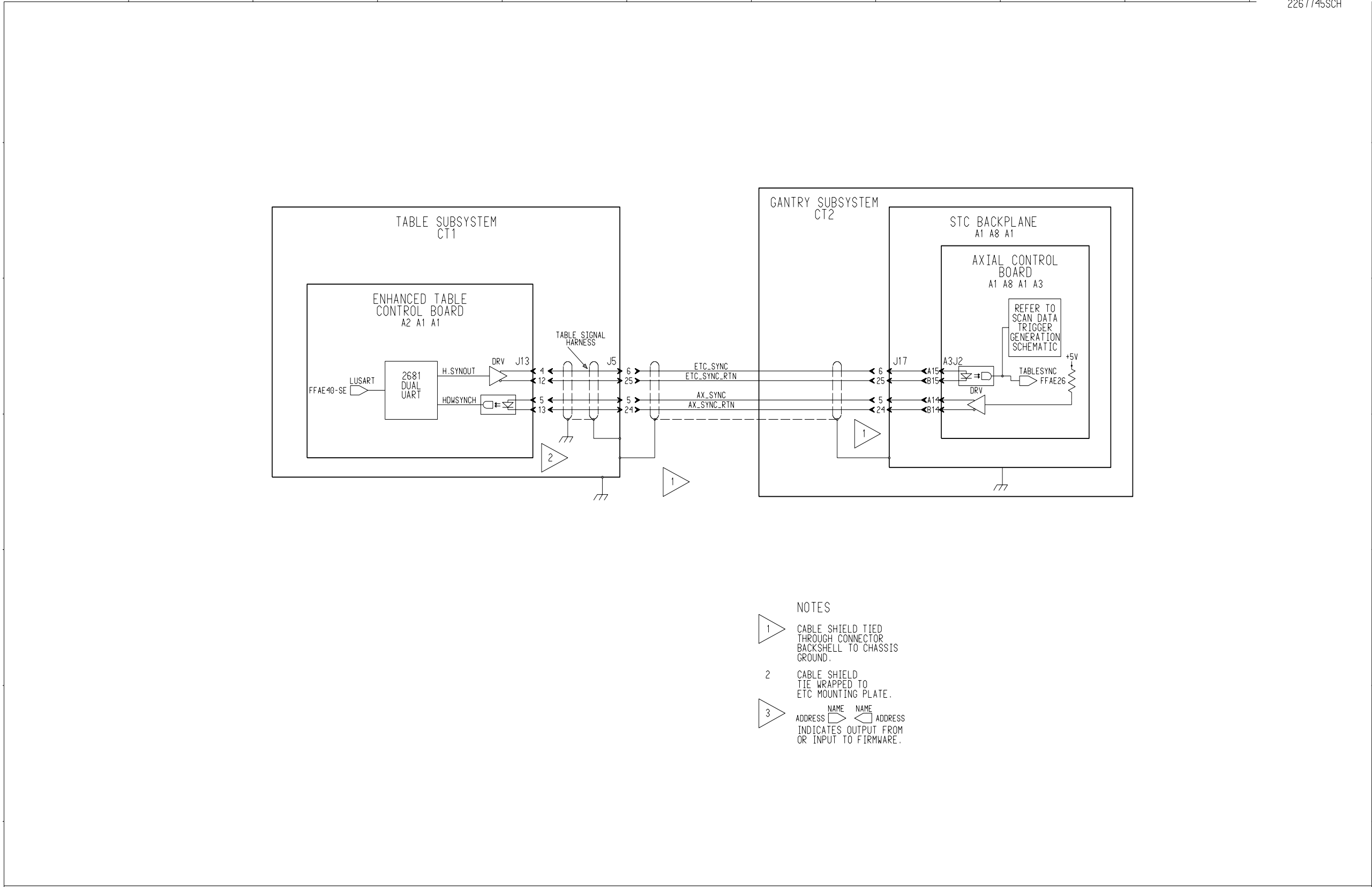


Figure 4-2 2267745 Sheet 1, Rev. 0

Section 3.0 — Table Elevation Drive - 2267754

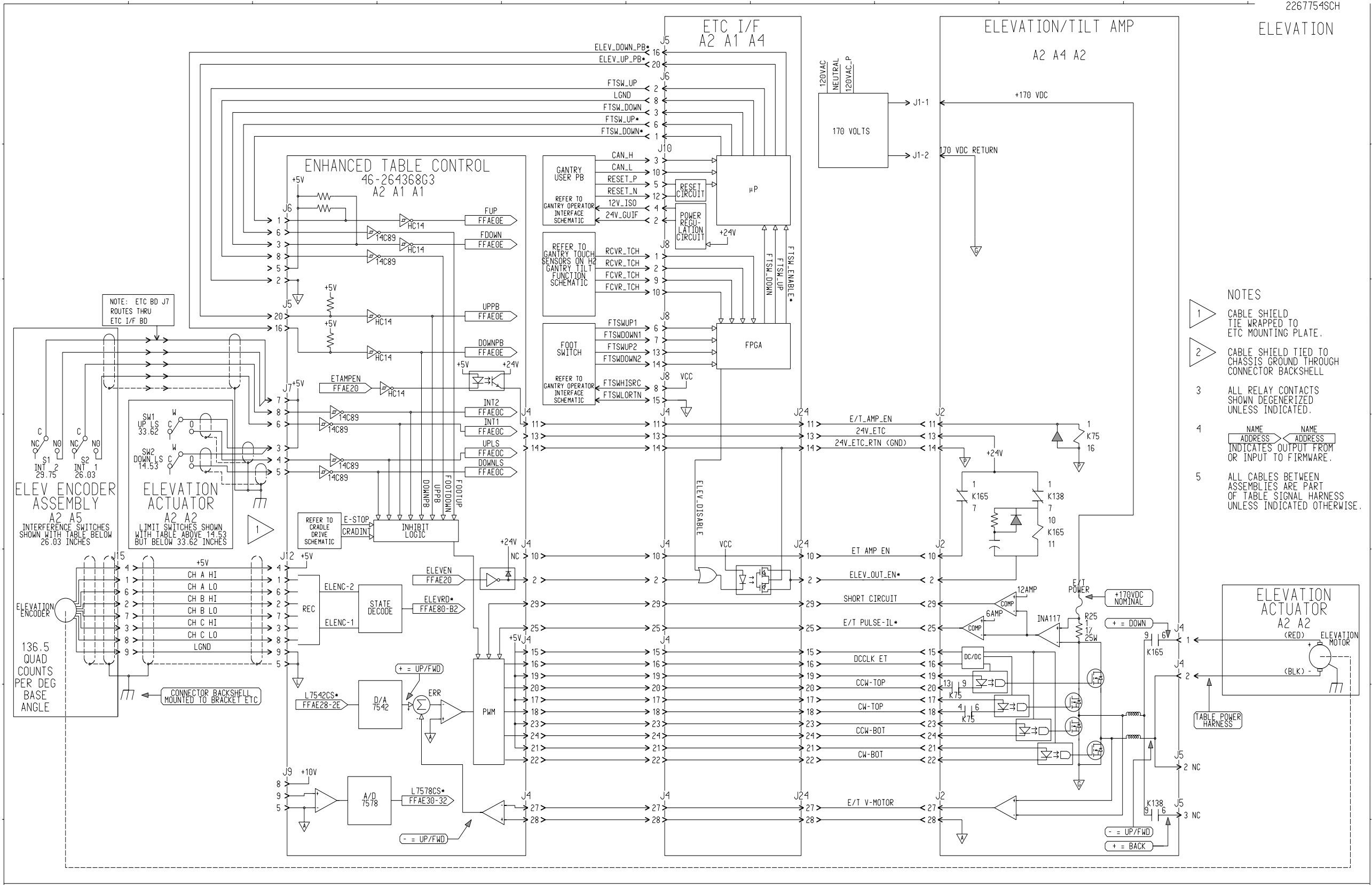

2267756SCH

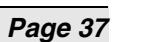

Section 5.0 — Table Cradle Drive (Vertical & Longitudinal Drives) - 2123540

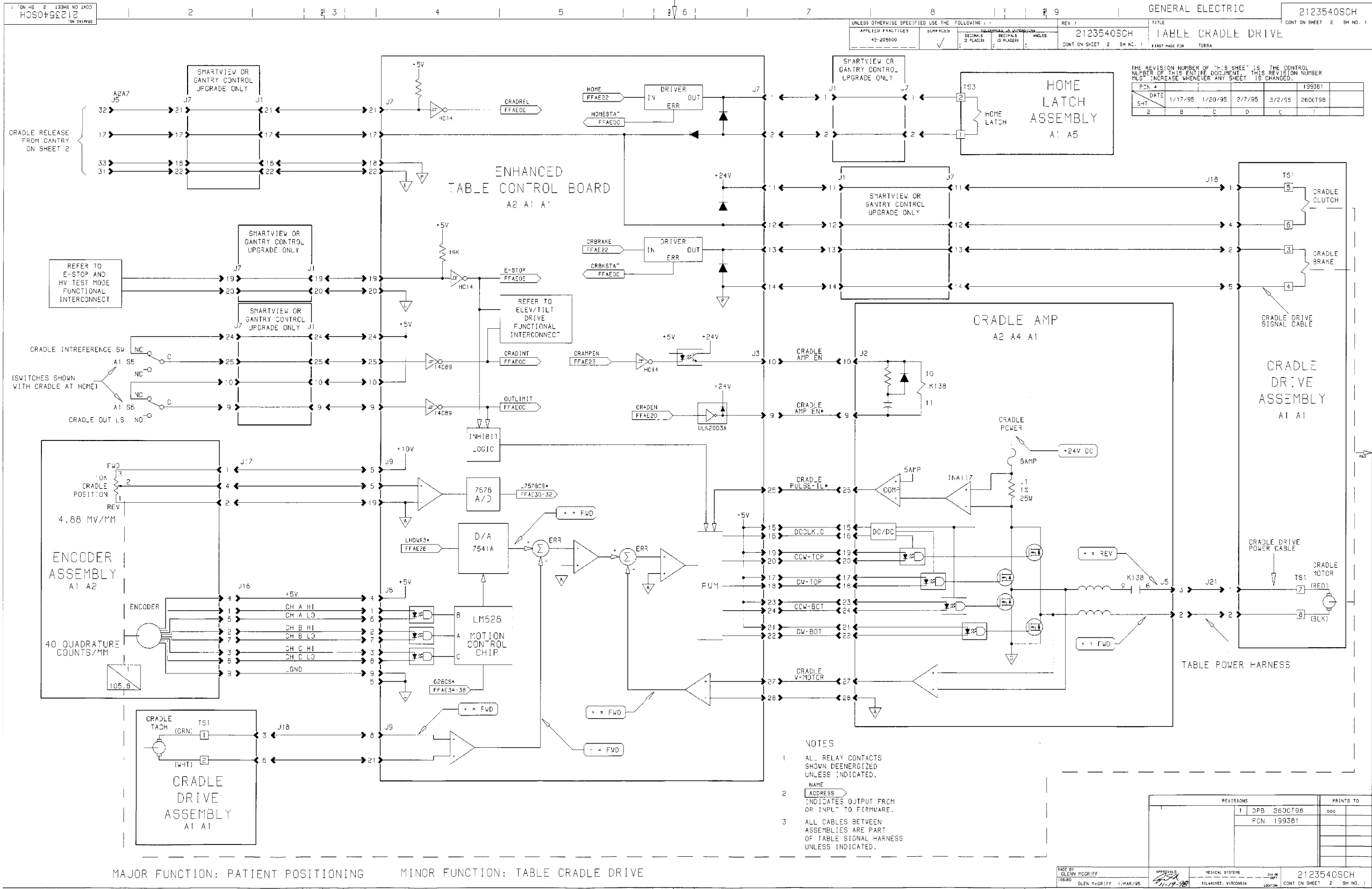


Figure 4-5 2123540 Sheet 1, Rev 1

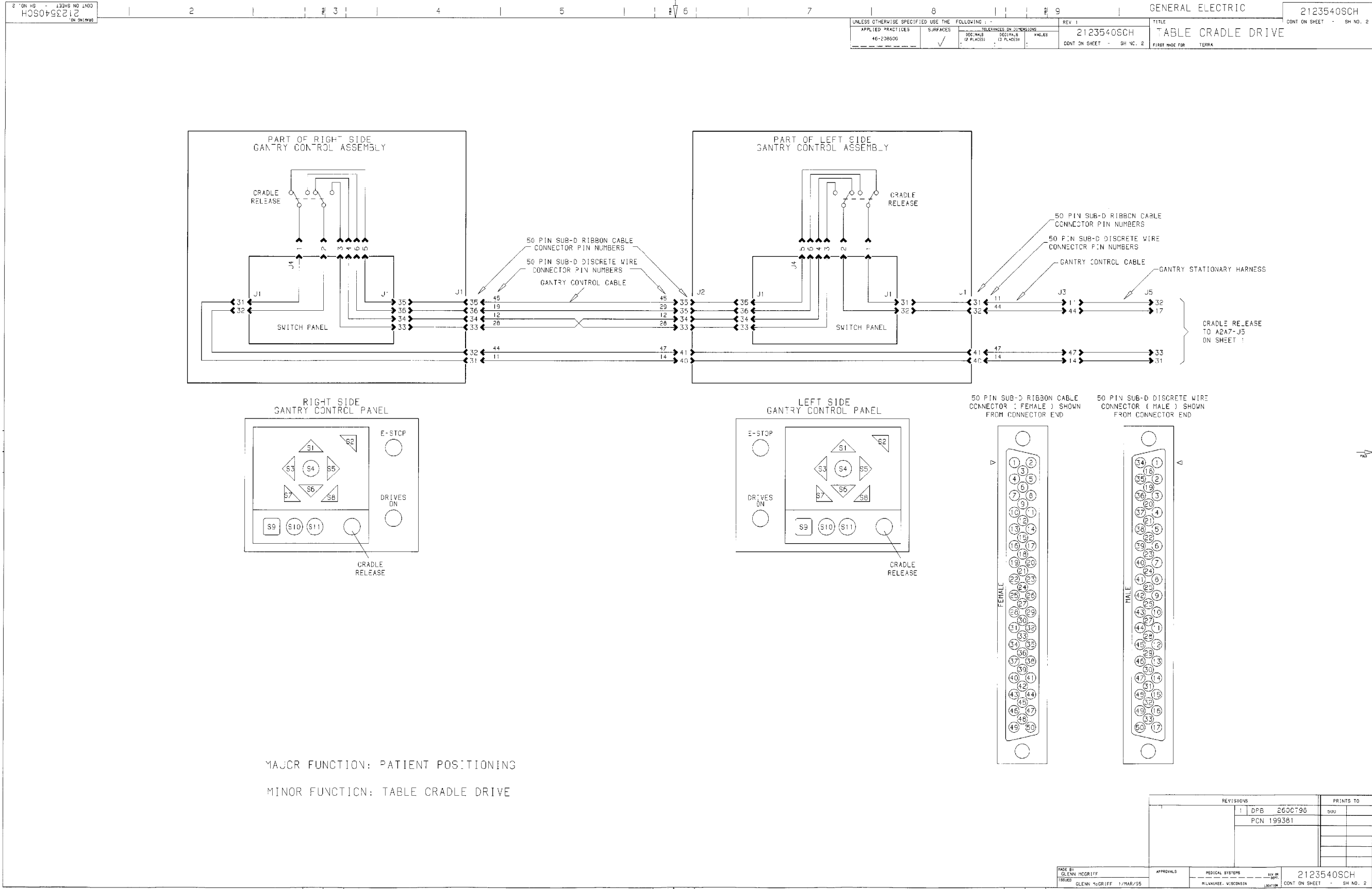



Figure 4-6 2123540 Sheet 2, Rev 2

NOTES:

Chapter 5

Axial Control

226 / 15 / SCH

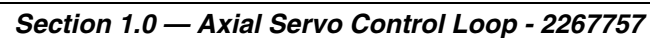


NOTES:

Chapter 6

Operator I/O

Section 1.0 — Gantry Operator Interface - 2267760

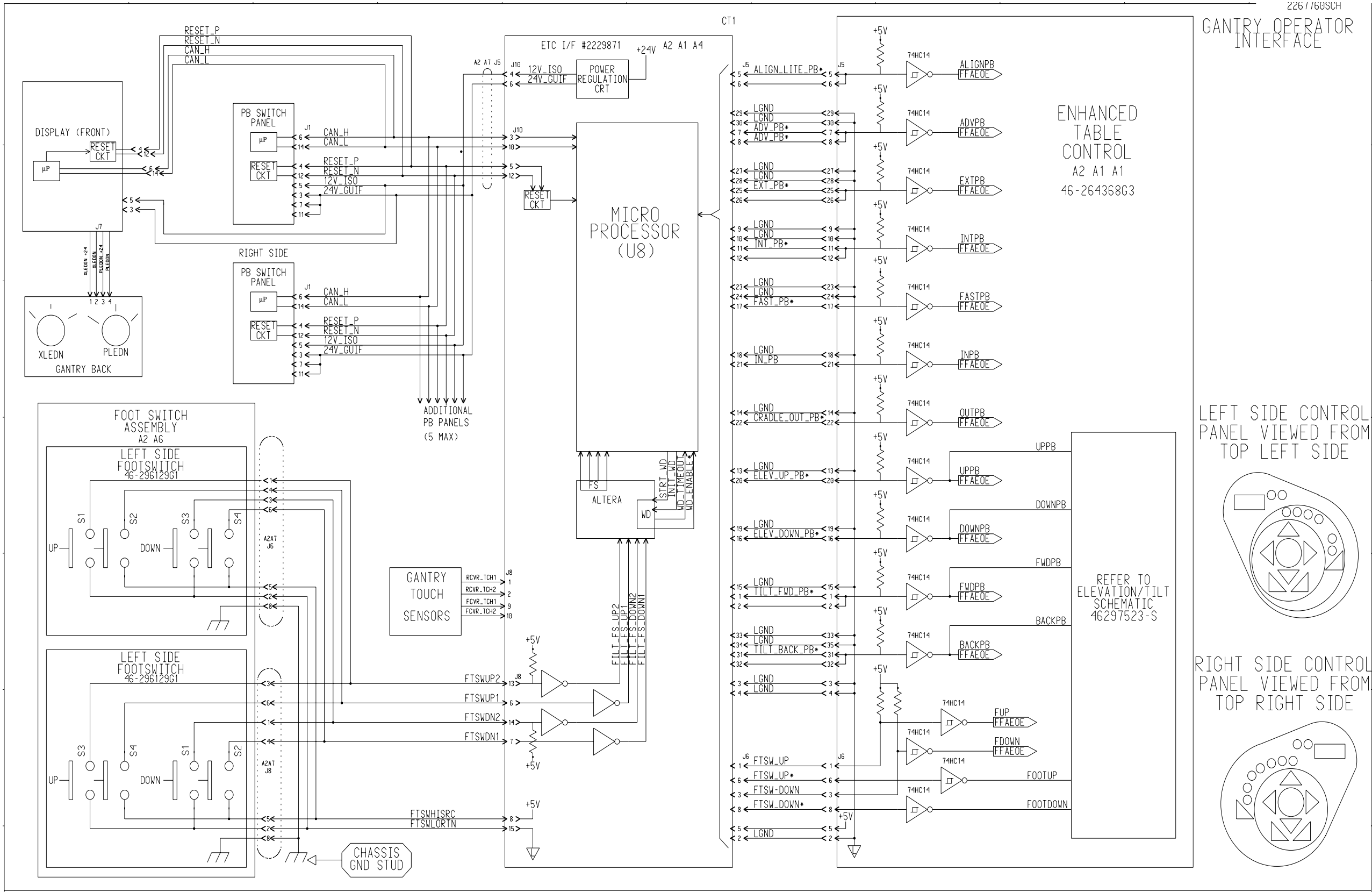

Section 2.0 — Site X-Ray ON Light - 2267761

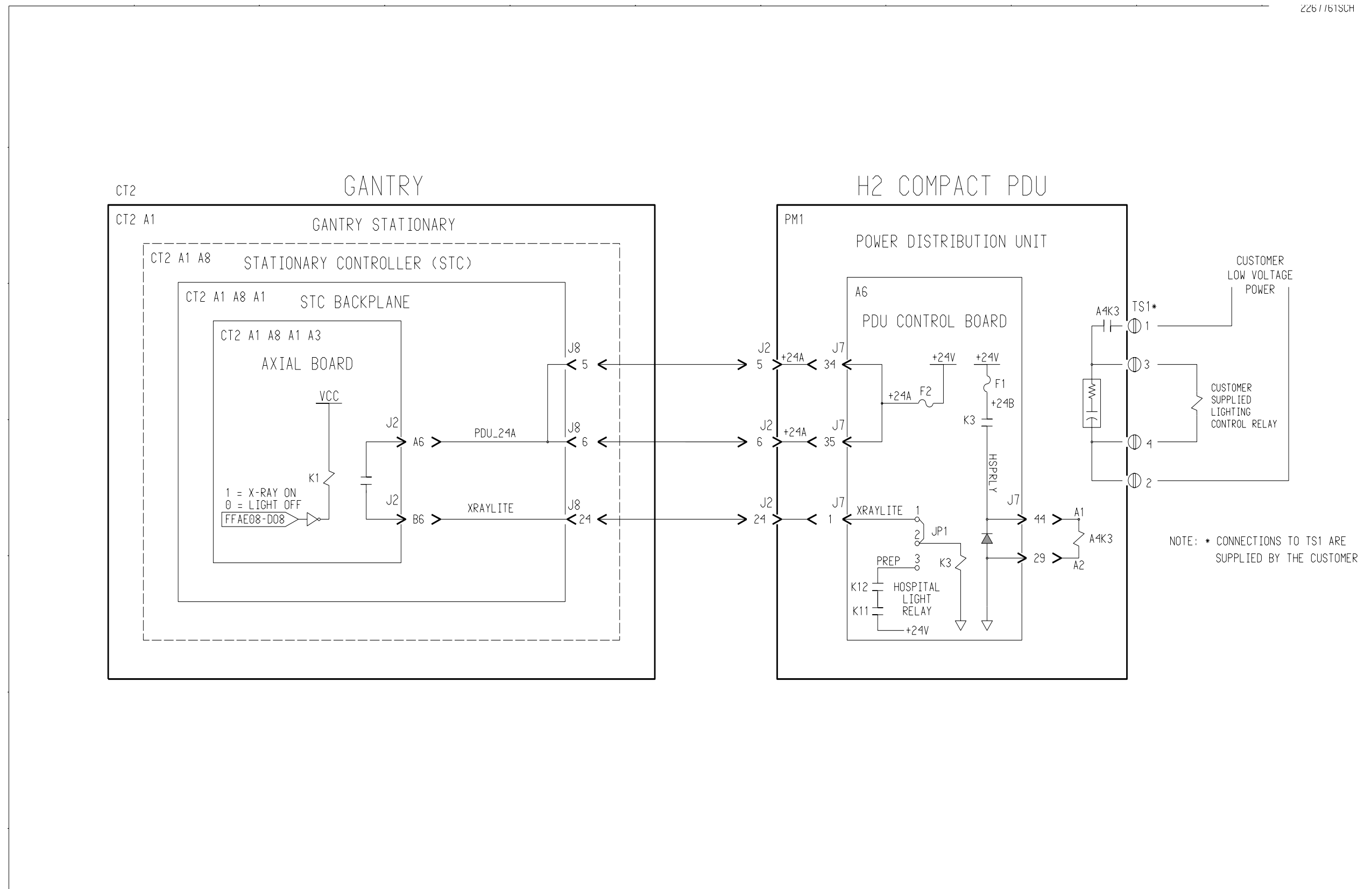


Figure 6-2 2267761 Sheet 1, Rev. 0

Section 3.0 — Gantry Options - 2286501

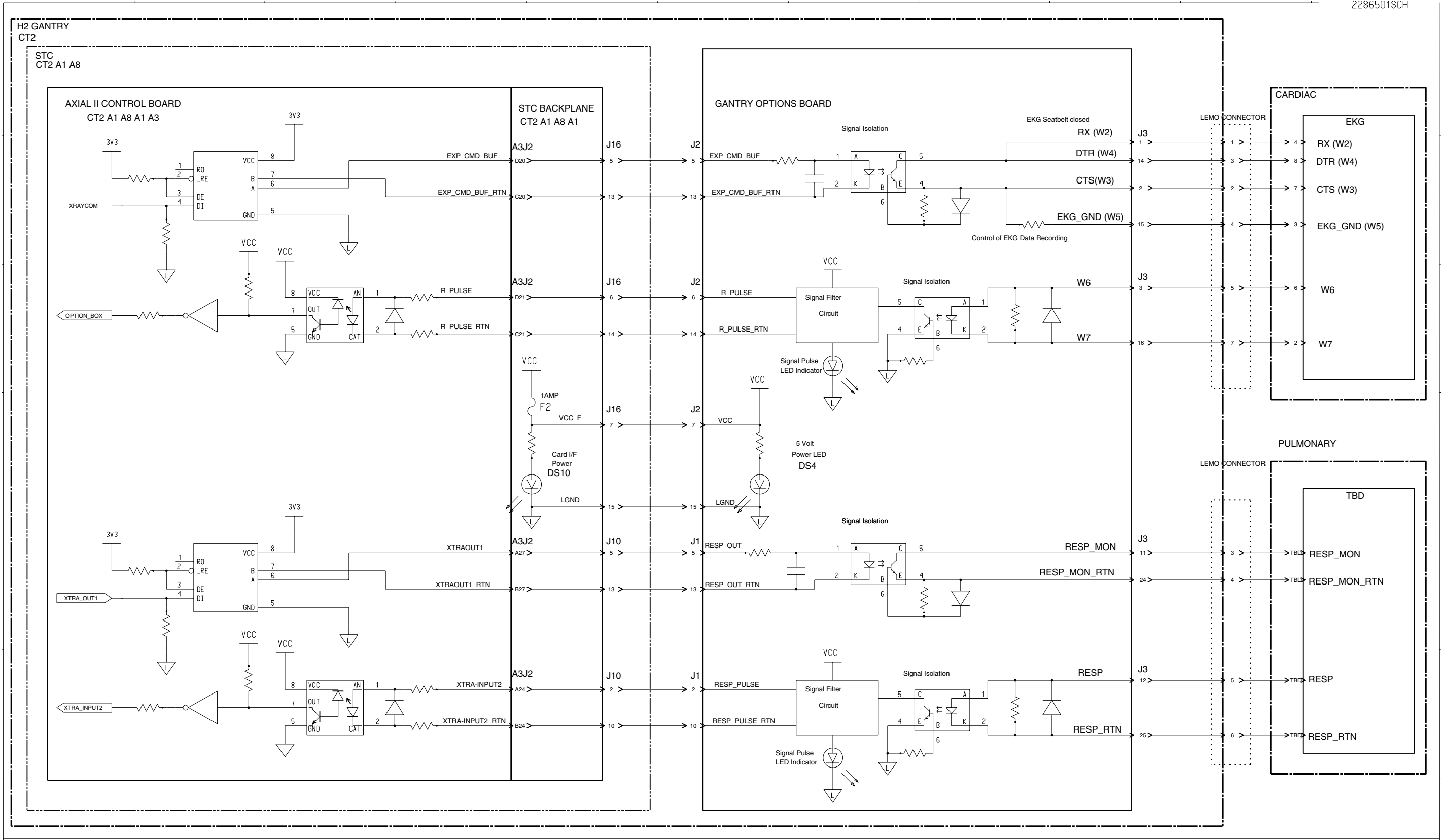


Figure 6-3 2286501 Sheet 1, Rev. 0

Section 4.0 — Smartview Controls - 2286504

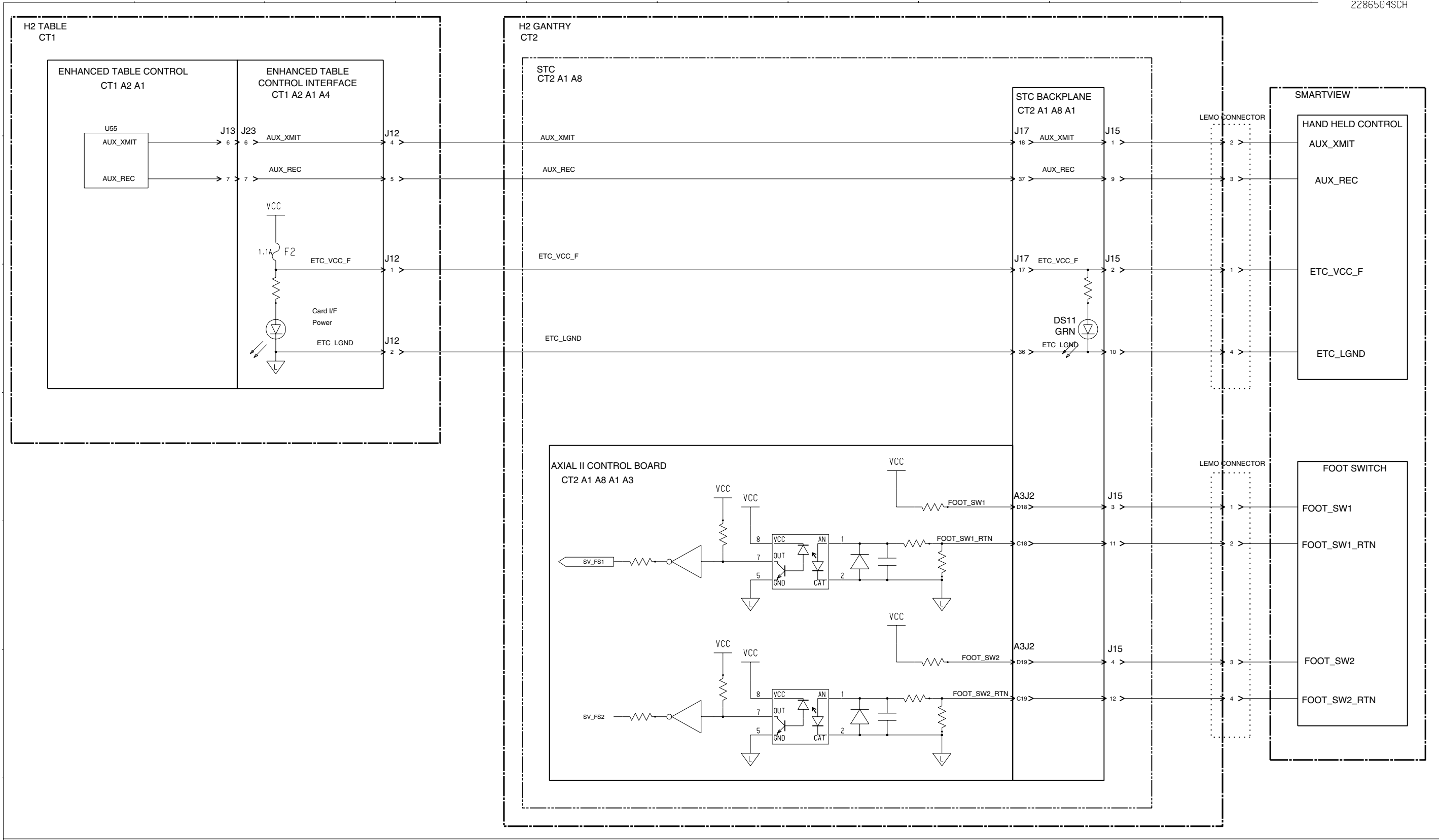


Figure 6-4 2286504 Sheet 1, Rev. 0

Section 5.0 — RX Tilt & Gantry Console Reset - 2267765

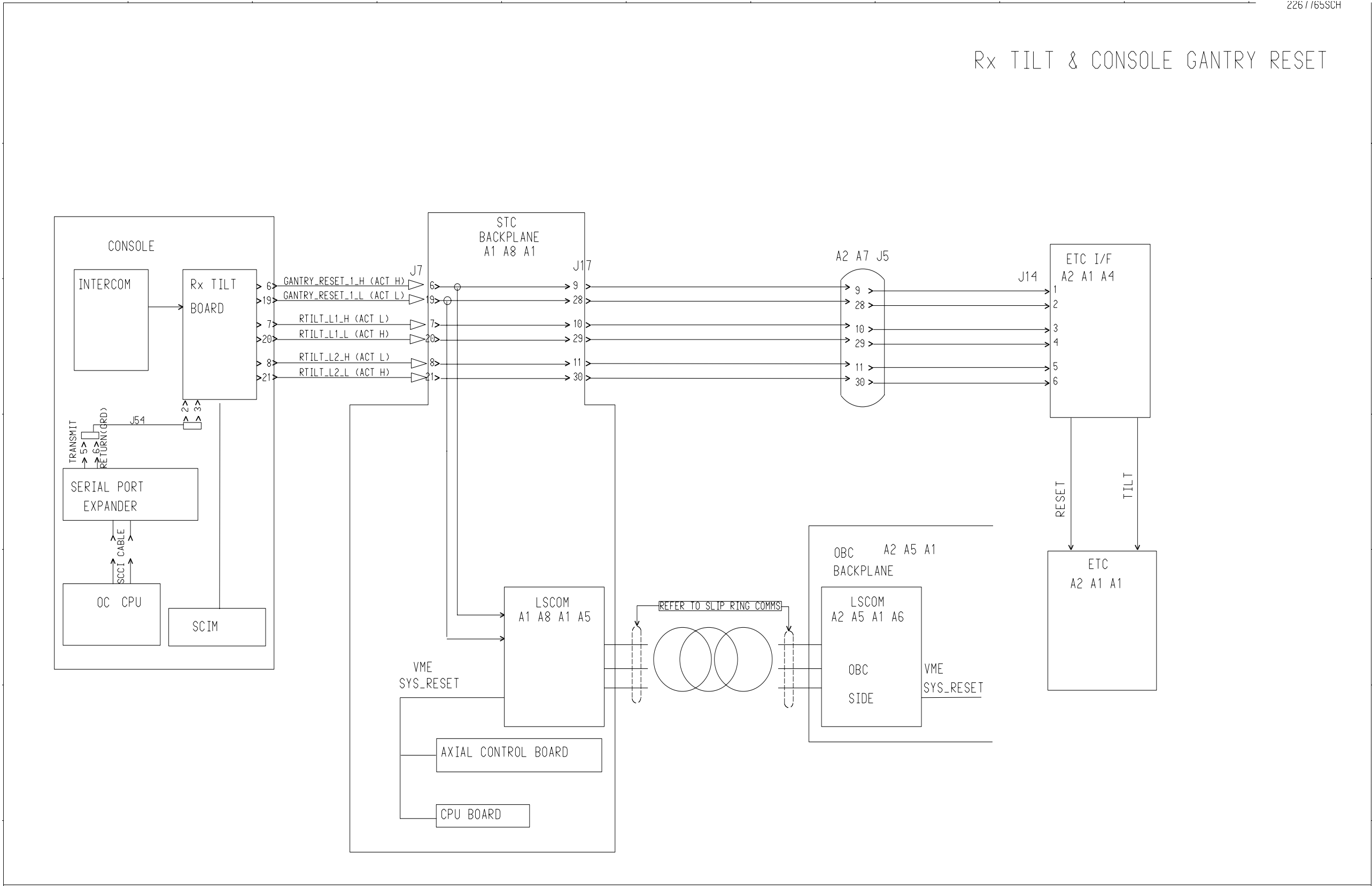


Figure 6-5 2267765 Sheet 1, Rev. 0

Chapter 7

System Monitoring

Section 1.0 — Mains Undervoltage - 2267731

2267731SCH

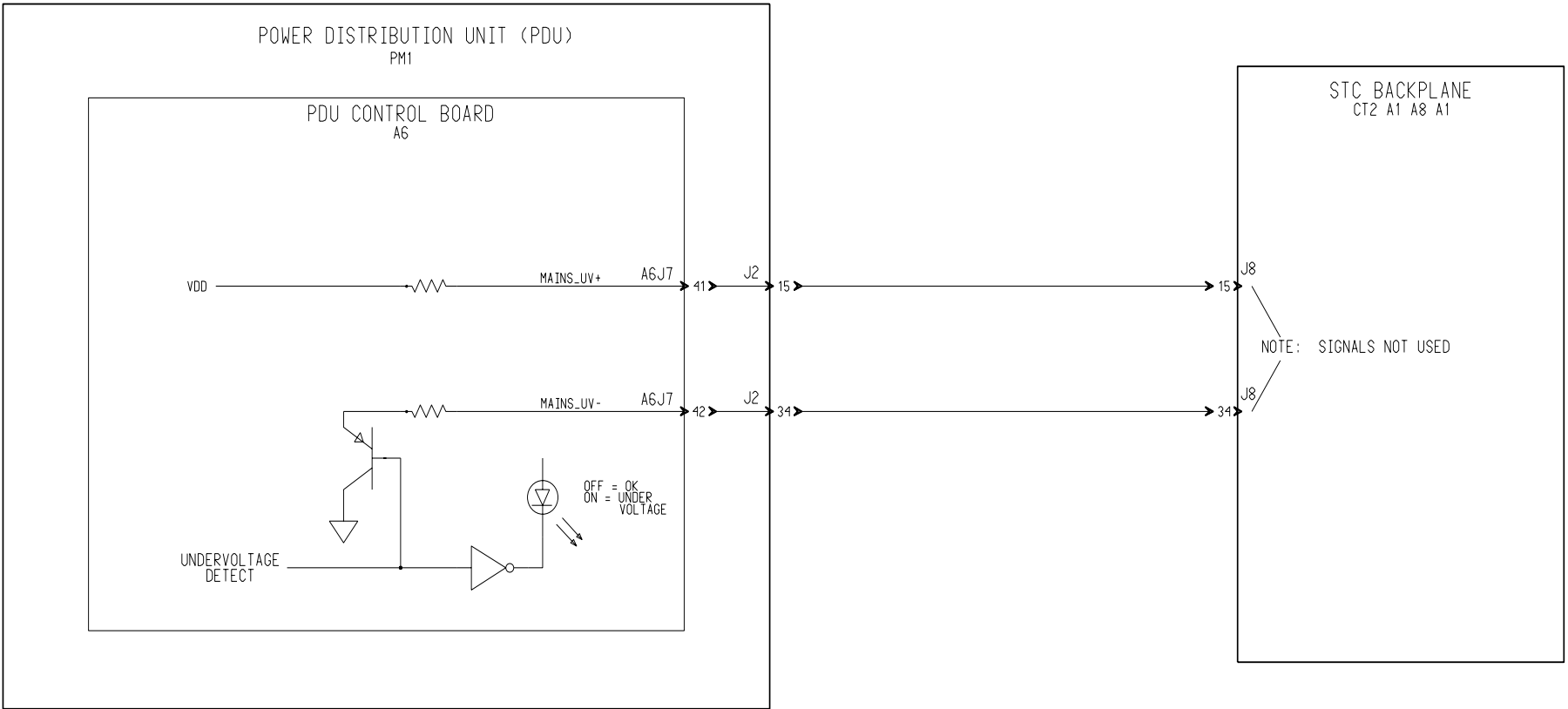


Figure 7-1 2267731 Sheet 1, Rev. 0

Section 2.0 — Gantry Temp. Sensors & Power Supply Monitors - 2267766

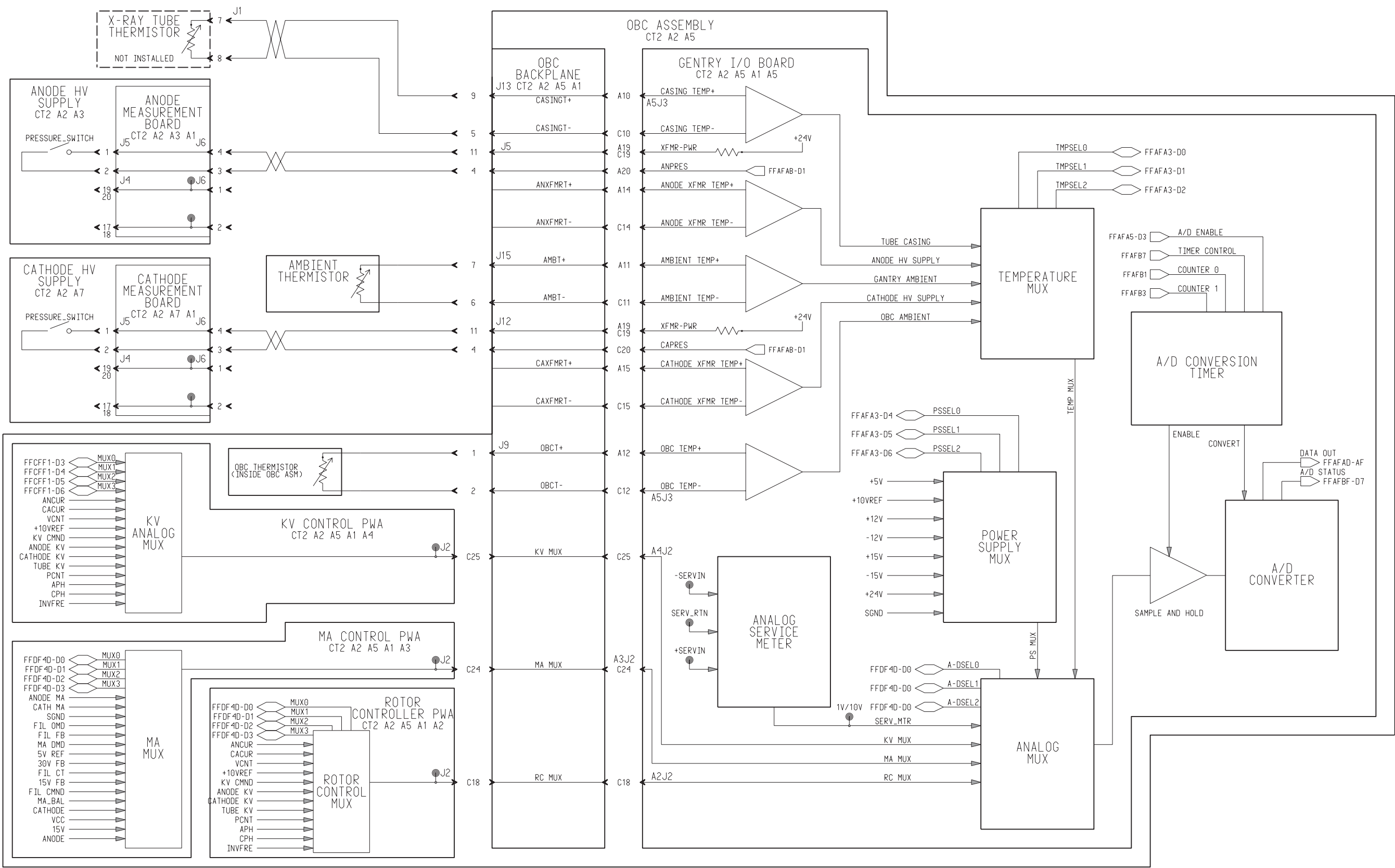


Figure 7-2 2267766 Sheet 1, Rev. 0

NOTES:

Chapter 8

System Power Control

Section 1.0 — HVDC Bus Control - 2267769

226 / 169SCH

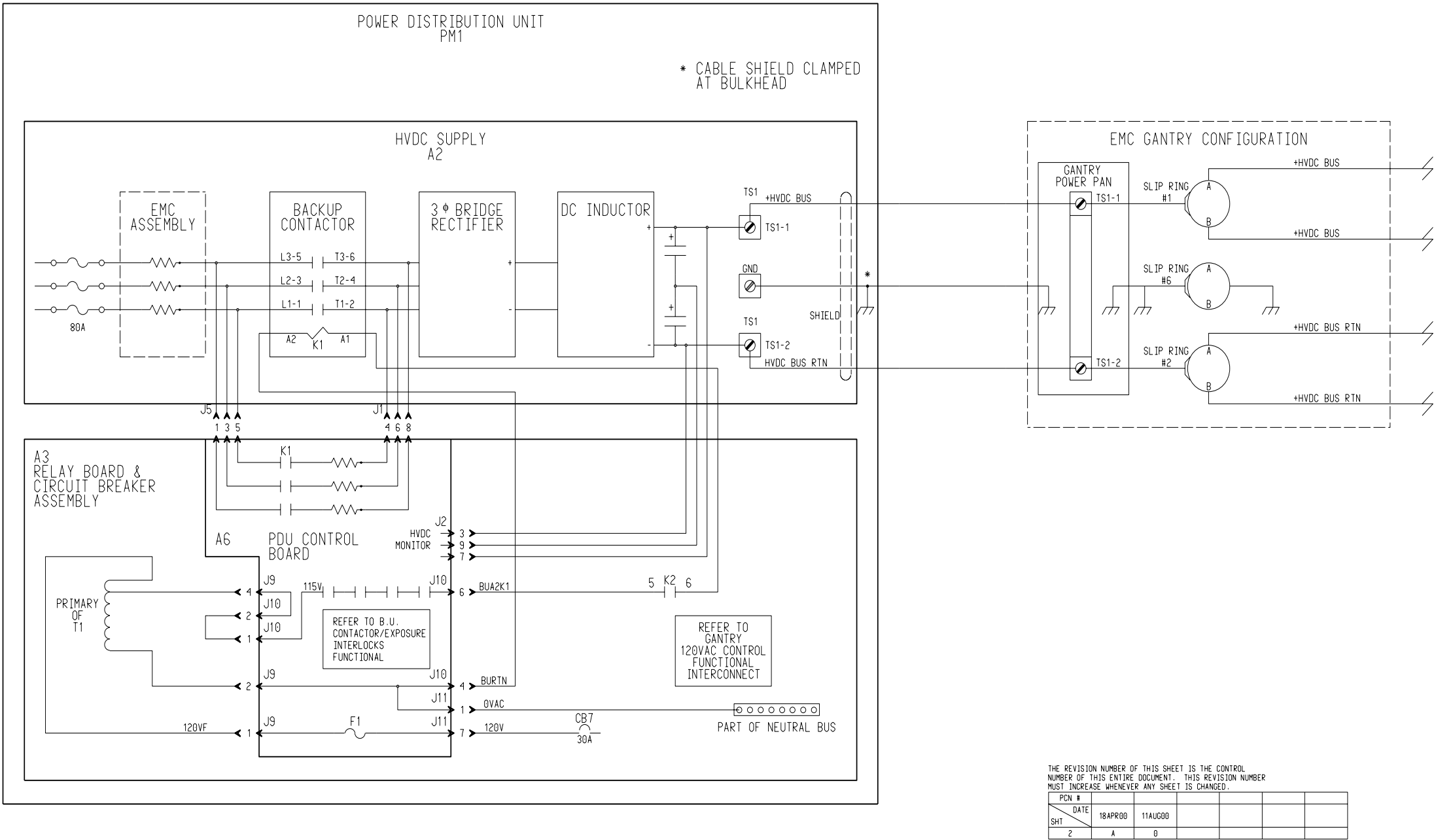


Figure 8-1 2267769 Sheet 1, Rev. 0

2267769SCH

SYSTEMS WITH PERFORMIX TUBE

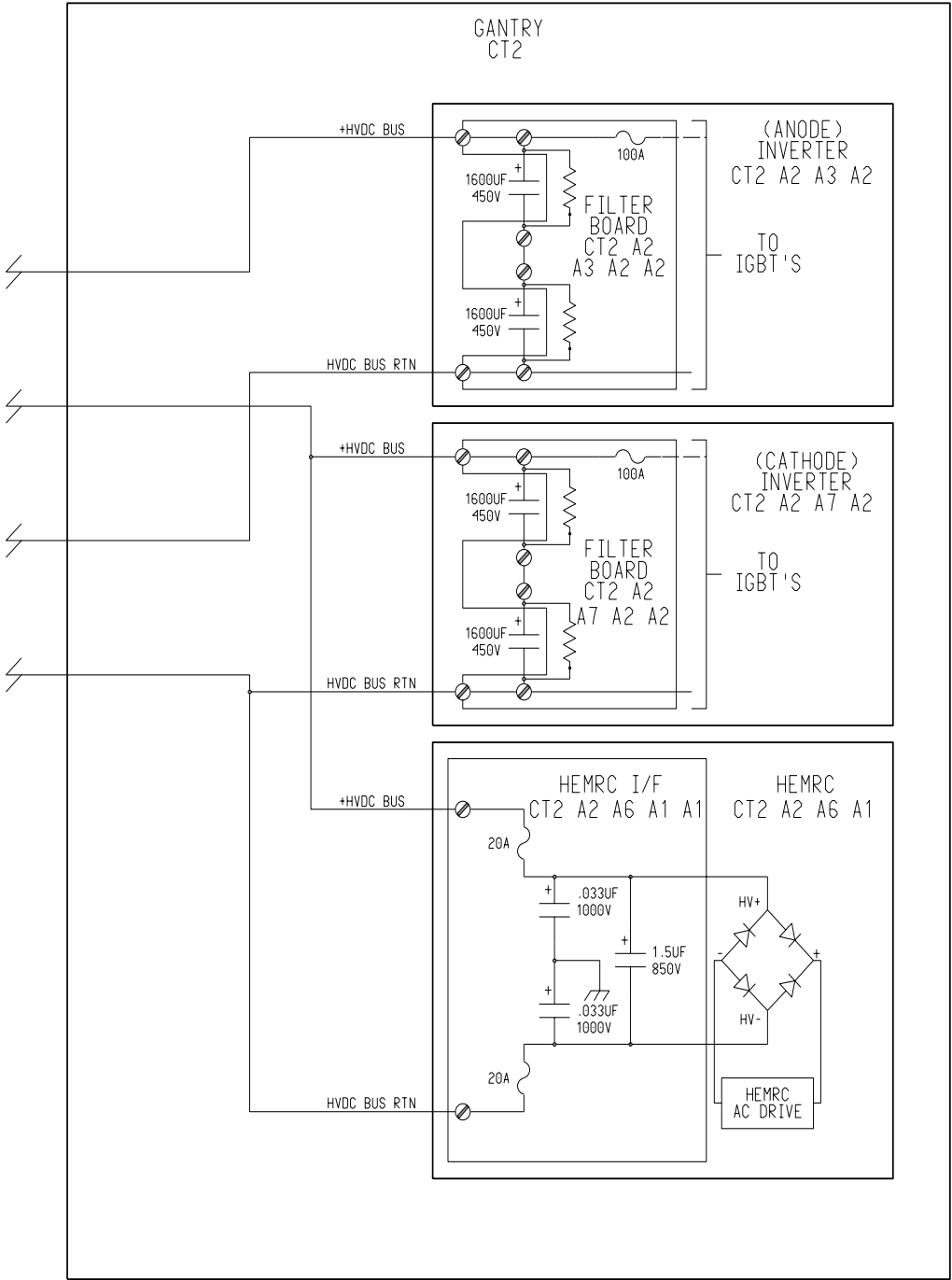


Figure 8-2 2267769 Sheet 2, Rev. 0

Section 2.0 — Gantry 120VAC Control - 2267770

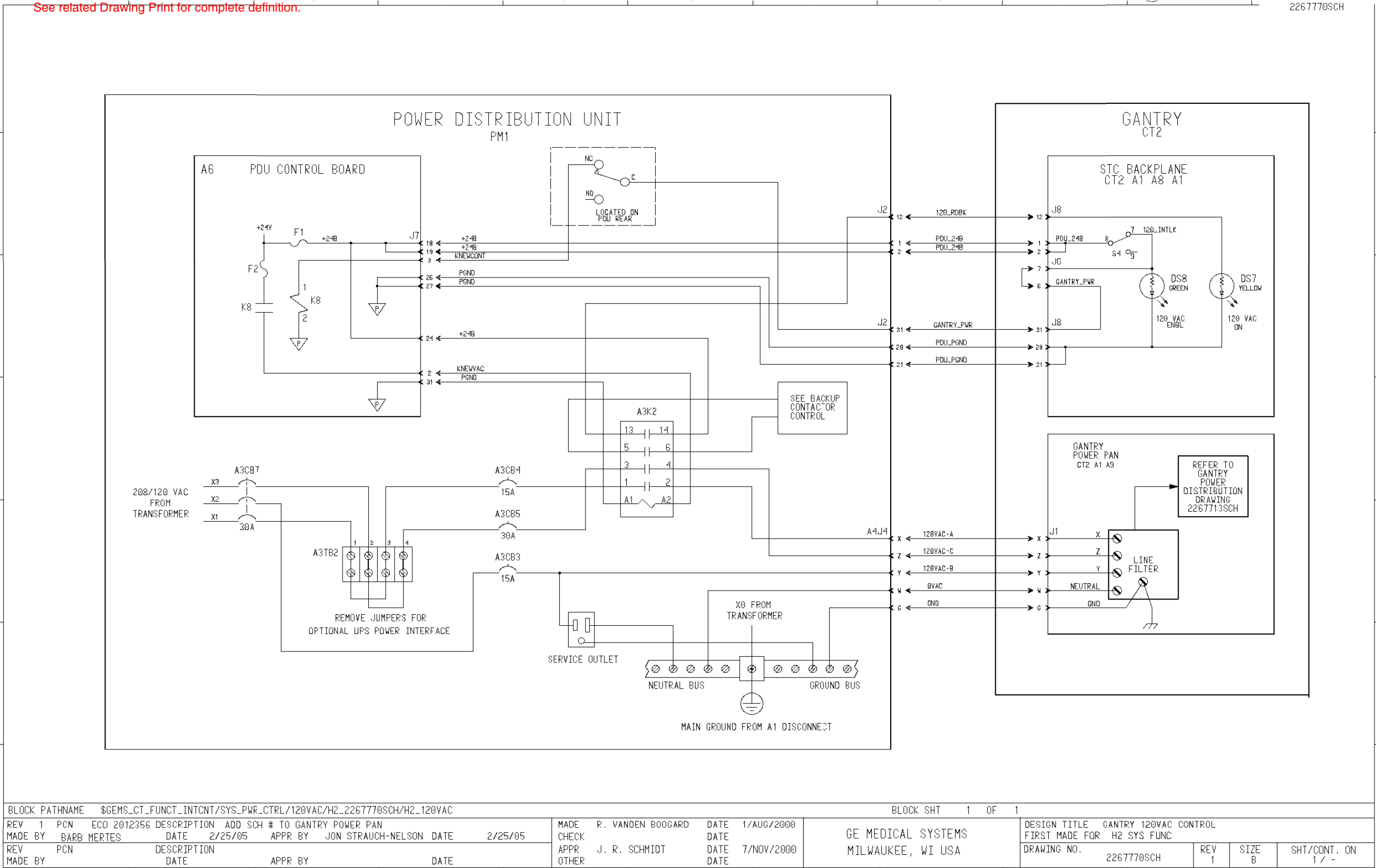


Figure 8-3 2267770 Sheet 1, Rev. 1

Section 3.0 — E-Stop Drives Reset - 2267773

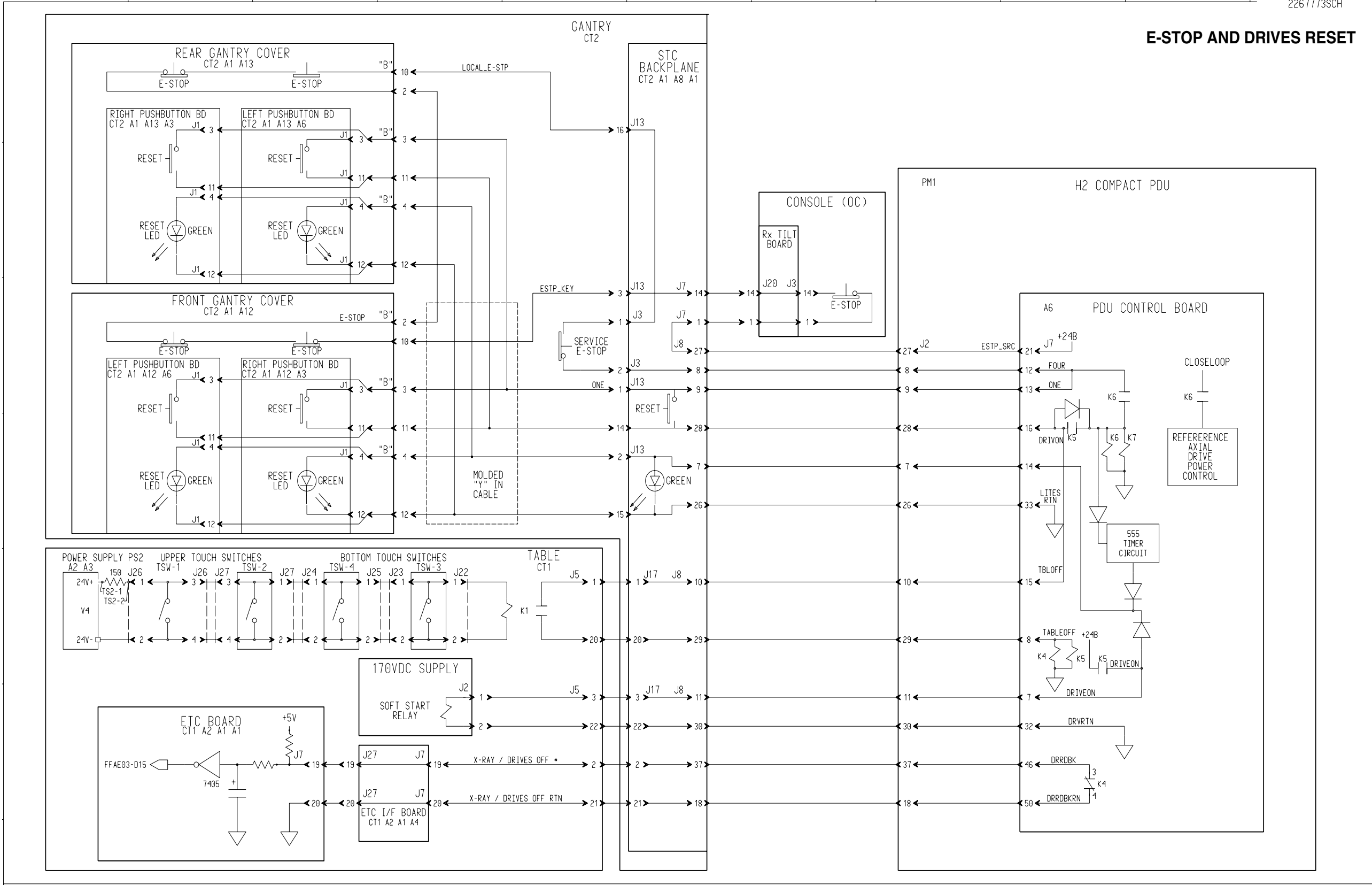


Figure 8-4 2267773 Sheet 1, Rev. 1

Section 4.0 — Gantry Stationary Power and Ground - 2267713

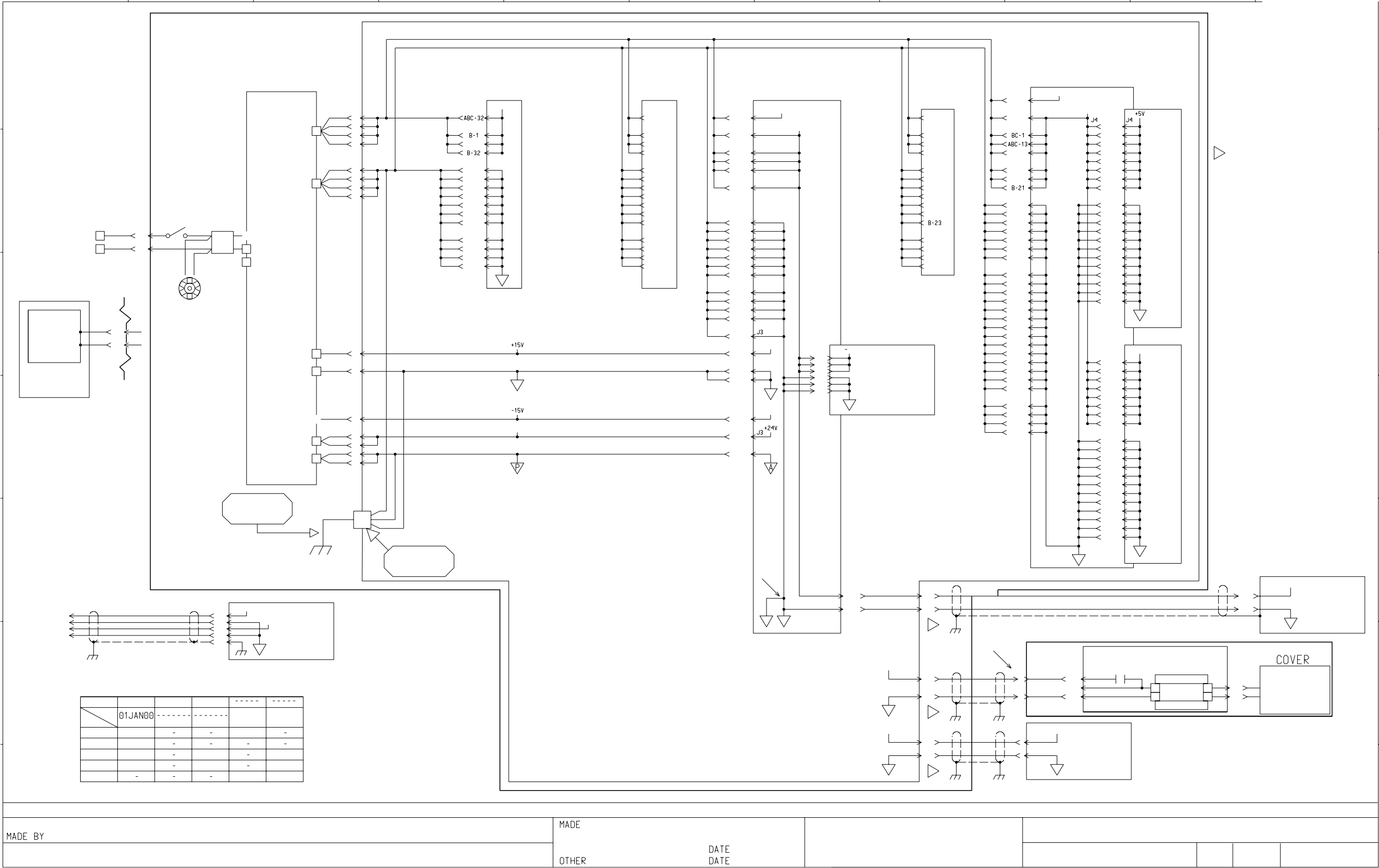


Figure 8-5 2267713 Sheet 1, Rev. 1

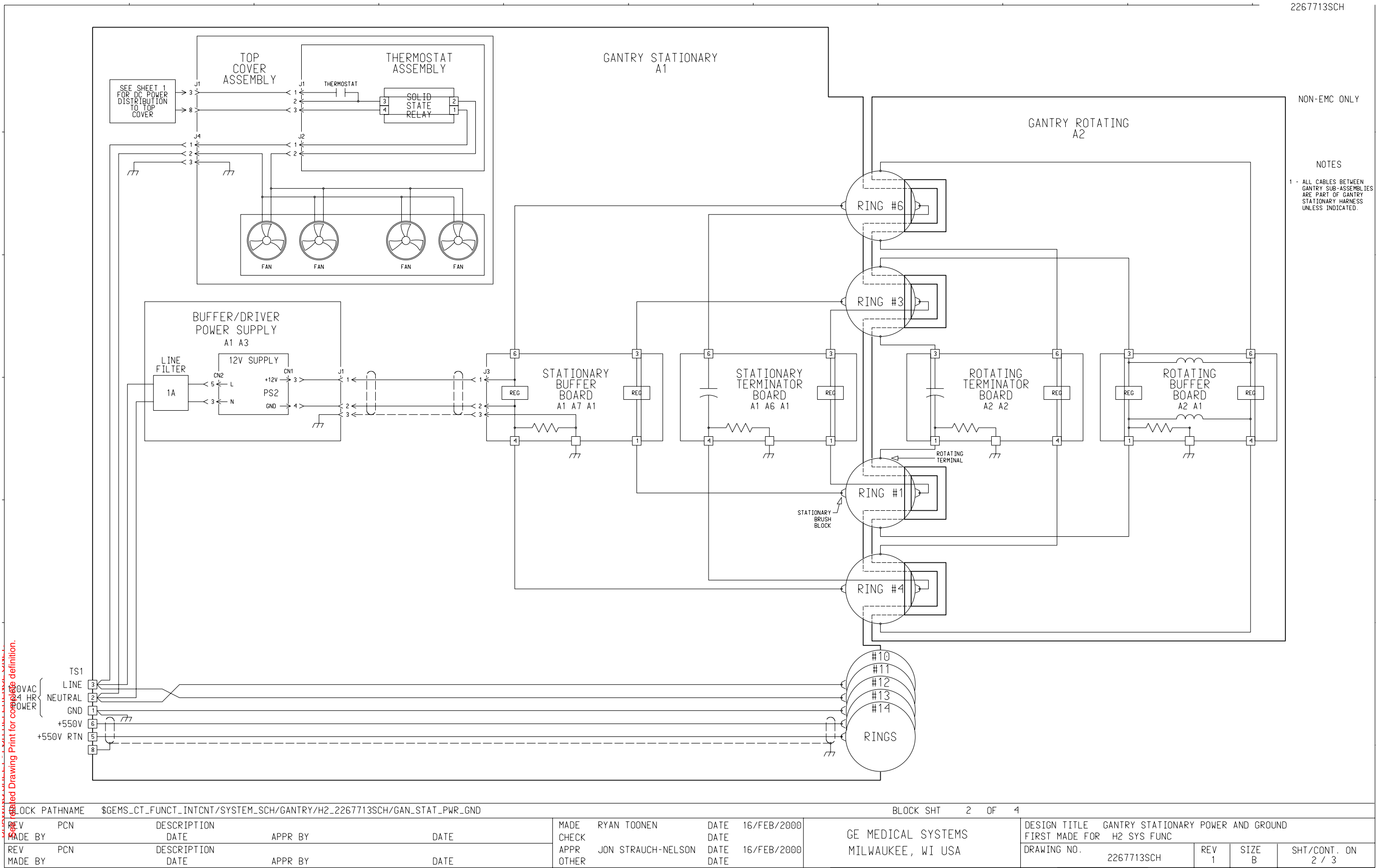


Figure 8-6 2267713 Sheet 2, Rev. 1

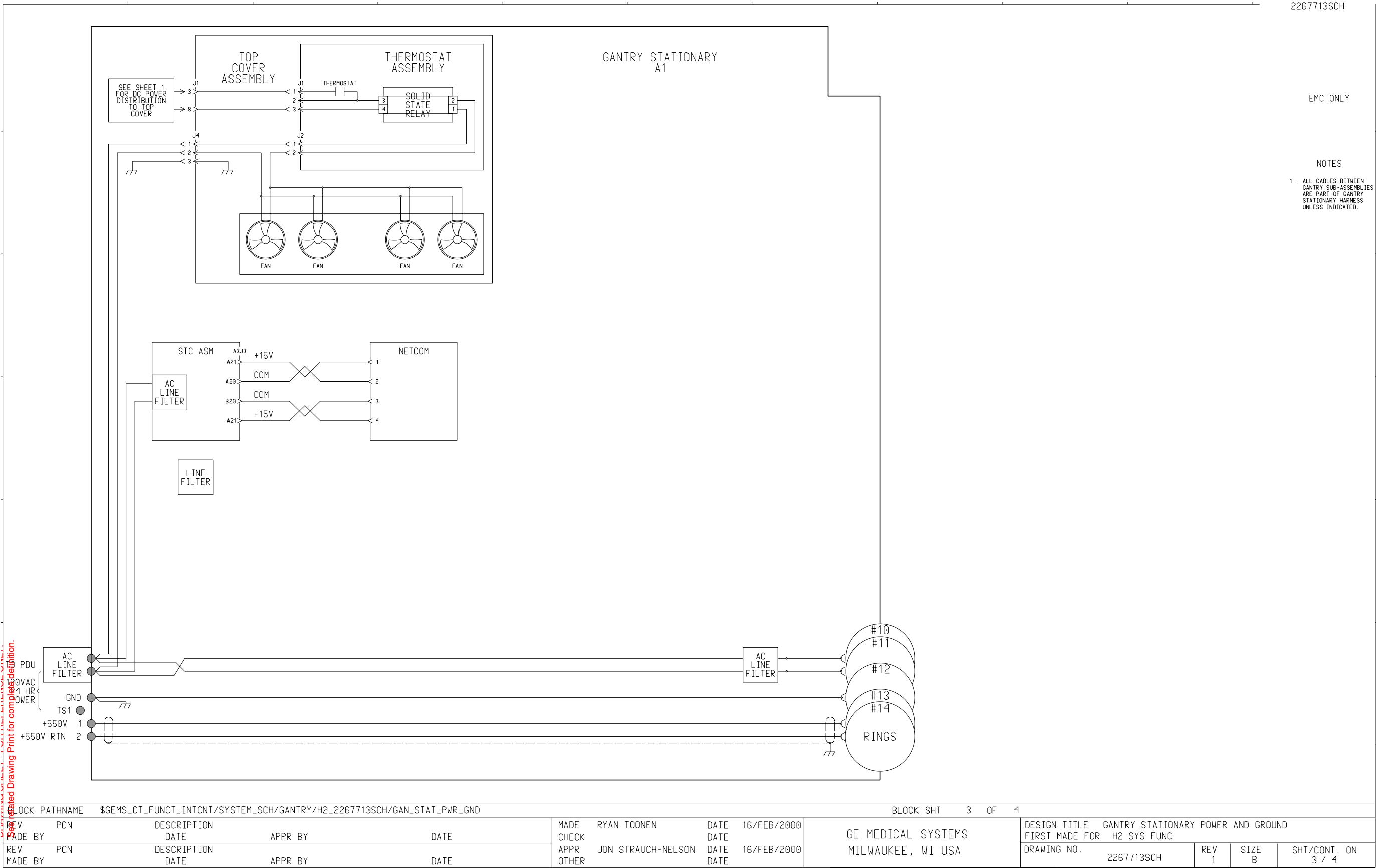


Figure 8-7 2267713 Sheet 3, Rev. 1




Page 63

NOTES:

Chapter 9

Image Generation

LightSpeed Plus Console (Signal Paths)

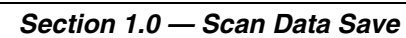

Section 2.0 — Data Restore

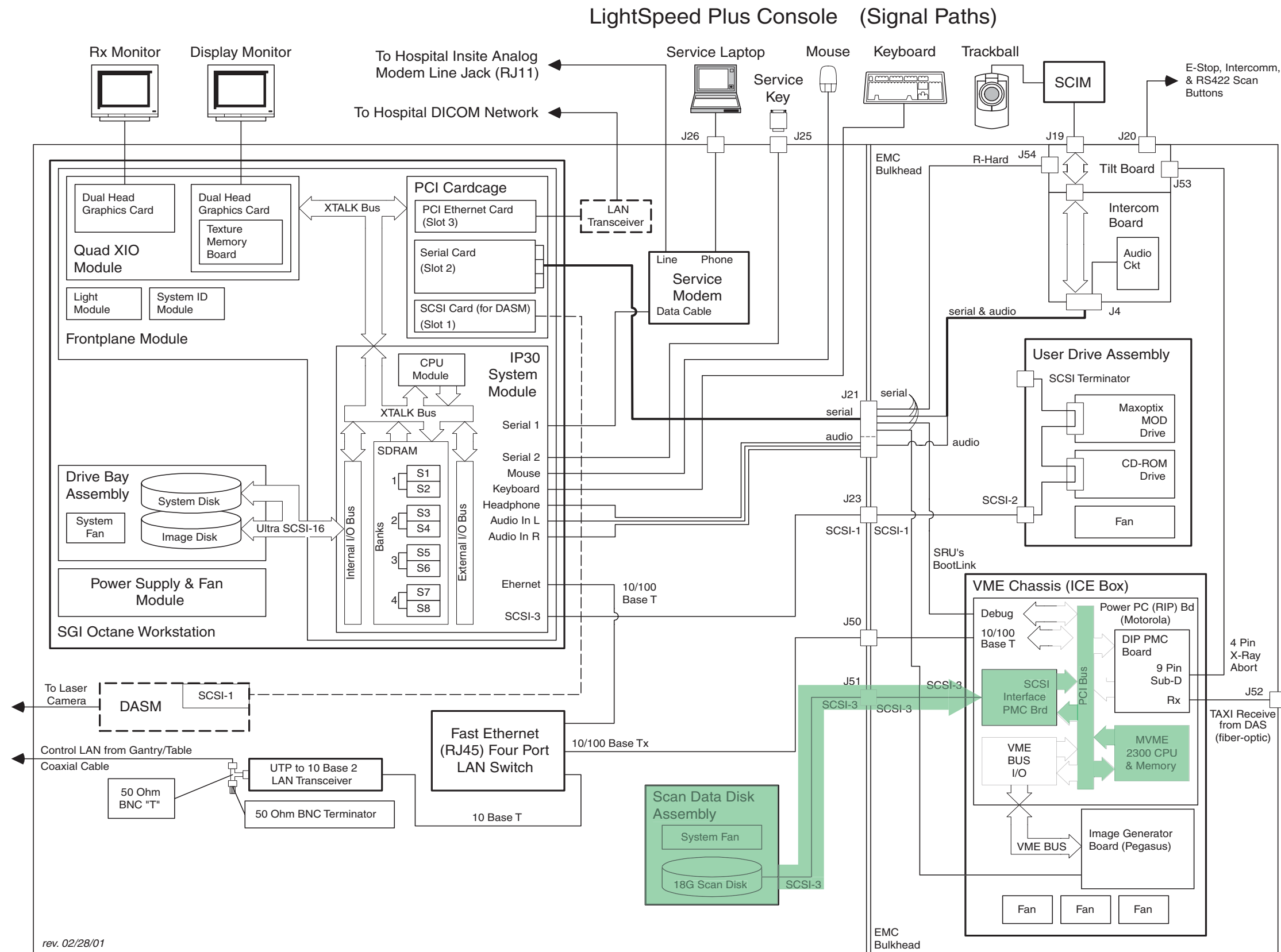


Figure 9-2 Data Restore

Section 3.0 — Pre-Processing, Filtered Back Projection & Post Processing

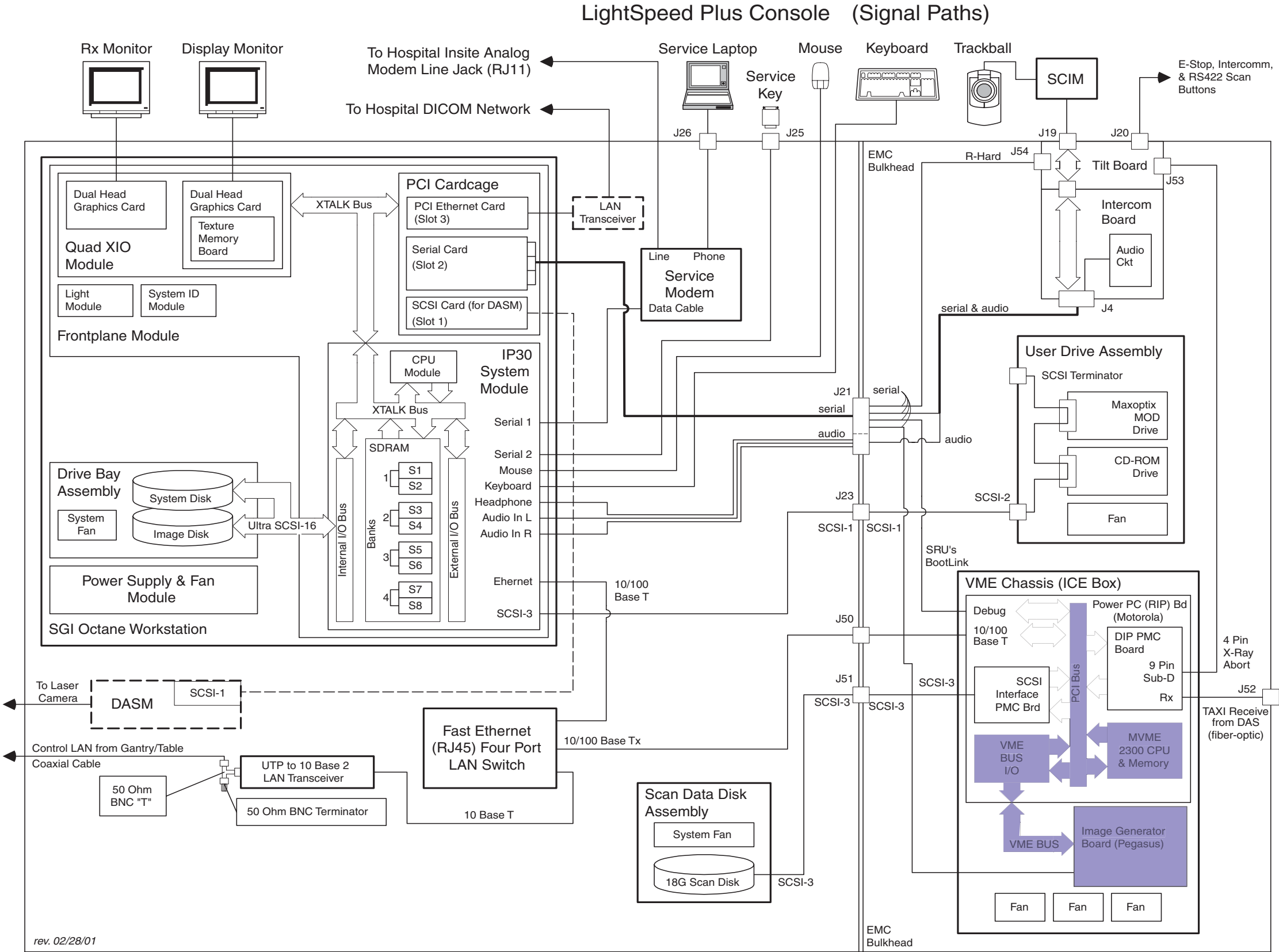


Figure 9-3 Pre-Processing, Filtered Back Projection & Post Processing

Section 4.0 — Image Transfer

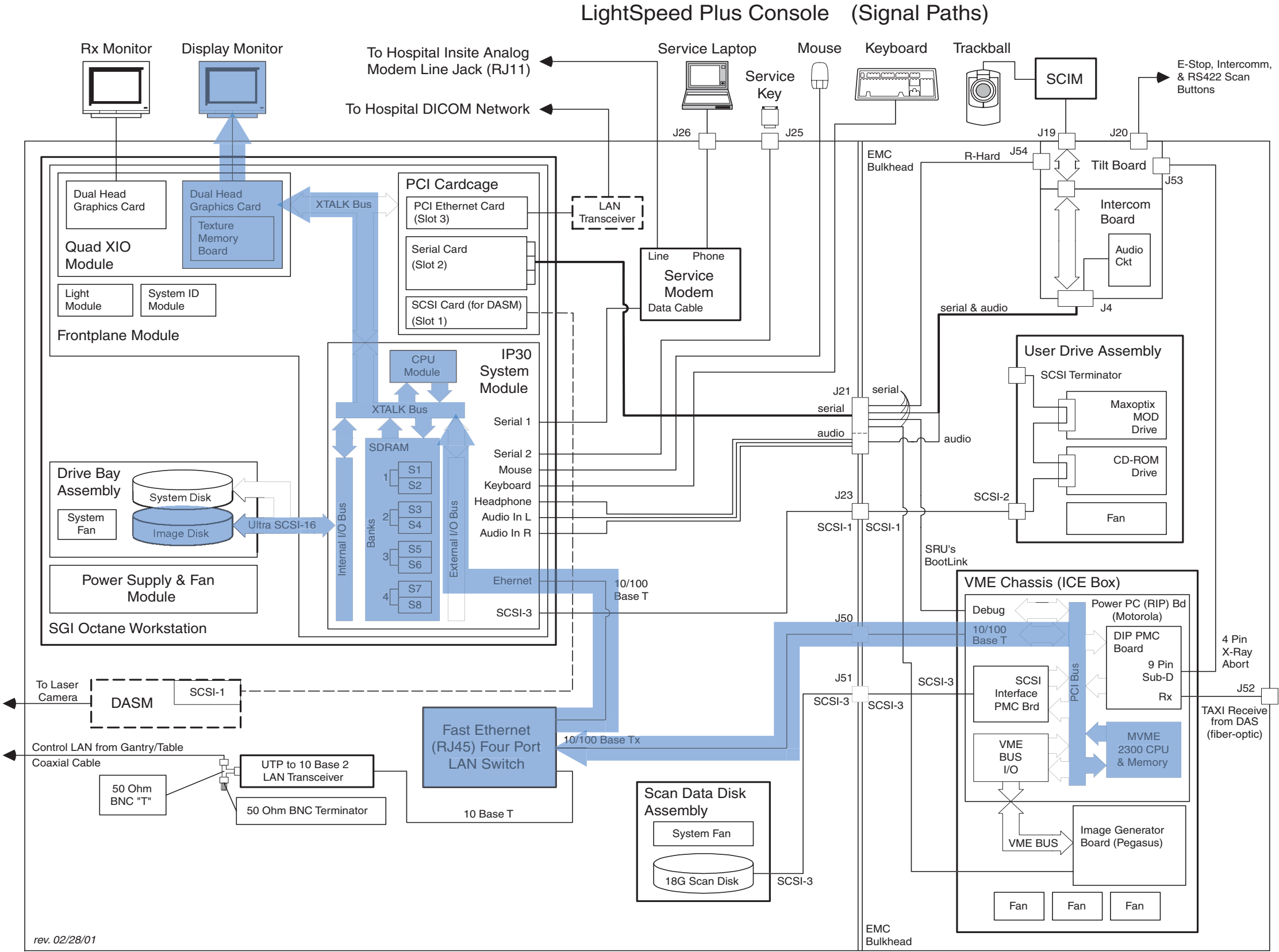


Figure 9-4 Image Transfer

NOTES:



GE MEDICAL SYSTEMS

GE MEDICAL SYSTEMS-AMERICAS: FAX 262.312.7434
3000 N. GRANDVIEW BLVD., WAUKESHA, WI 53188 U.S.A.
GE MEDICAL SYSTEMS-EUROPE: FAX 33.1.40.93.33.33
PARIS, FRANCE
GE MEDICAL SYSTEMS-ASIA: FAX 65.291.7006
SINGAPORE