



GE Medical System

Technical Publications

2321165-100
Revision 0

CT Mobile Startup/Shutdown Operator Manual

Do not duplicate

Copyright© 2001, General Electric Co.



LightSpeed Mobile Startup/Shutdown

Feature Description:

The LightSpeed Plus Mobile Startup and Shutdown procedure provides capability for the user to change site name as part of the start up procedure and ready the system for transport as part of the shut down procedure. The LightSpeed Plus Mobile Startup and Shutdown procedures will need to be followed each time the mobile is moved to a new location.

Start Up Procedure:

1. Complete the van manufacture's instructions for powering up and leveling the CT mobile van on the mobile pad.
2. Insure the system power to the CT scanner is turned on. Two green lights on the main power panel should be illuminated. If not, follow the "Start-Up" instructions posted on the Power Distribution panel.
3. Press the [Power On] switch on the console to begin the start up process.
4. In a few minutes, the CT Mobile Site Setup screen will be displayed to provide selection of site information.

Site Names:

- If there is only one site listed on the Mobile site list, the system will automatically complete the start up procedure.
- If there are multiple site names listed in the Mobile site list, the Mobile site list pop-up will be displayed.

Mobile Site Setup Screen

Mobile Site Setup

CT Mobile Site Setup

Site Information:

Site Name:

SuiteID:

Hostname:

IP Address:

Netmask:

Broadcast:

Default Gateway:

Hospital A
Hospital B
Hospital C
Hospital F

Add Change Remove Search Select Done

Select desired site name from the site list and click on [Select].

Or

- Select [Search] to find a site by name or hostname.

Mobile Site Setup Screen

The screenshot shows the 'Mobile Site Setup' window. The main form is titled 'CT Mobile Site Setup' and contains a 'Site Information:' section with the following fields: Site Name, SuiteID, Hostname, IP Address, Netmask, Broadcast, and Default Gateway. To the right of these fields is a list of registered sites: Hospital A, Hospital B, Hospital C, and Hospital F. Below the Site Information section are buttons for 'Add' and 'Change'. A 'Search Mobile site' dialog is open in the foreground, showing a 'Search by' section with two input fields: 'Hospital' and 'Hostname'. The 'Search' and 'Done' buttons are at the bottom of the dialog.


Enter the site name in Hospital field or the network name in Hostname and select [Search] to search the list of registered sites for the designated site.

This screenshot shows the same 'Mobile Site Setup' window as the previous one, but the 'Search Mobile site' dialog is now in a different state. The 'Hospital' input field now contains the text 'Hospital B'. The 'Search' and 'Done' buttons remain visible at the bottom of the dialog.

Select [Done] to exit Search mode.

The system will display the Mobile Site confirmation pop-up.

Mobile Site Activation Screen

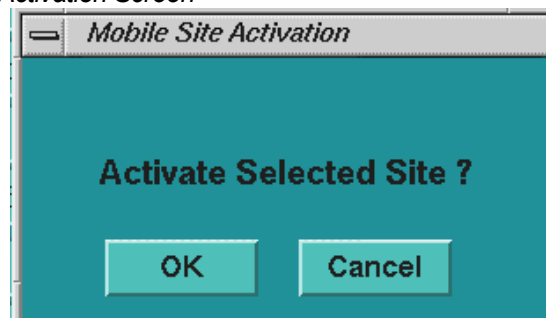

Select [OK] to accept the selected site and continue.


Select [Cancel] return to Mobile Site Setup screen.

- The site name and network information will be entered and the system start up will continue.
 - As the system starts up, there is a series of screens that appear. Please wait till you see the "Reference cradle and elevation drives." message in the current message area and the "OC Initializing" message disappears before continuing. These messages indicate the system is now ready for use.
5. Remove the Table Transport Pole from under the table. Store the pole in cabinet until needed to ready the system for transport.
 6. Press the [Table Down] button to lower the table to reference elevation. Lower the cradle until the system recognizes the table height. Press together the [Tilt Backward] and [Tilt Forward] buttons to bring the gantry to zero degree position. After completion of cradle reference, the system is ready.
 7. Perform Fast Cals every day. Perform a Tube Warm Up if Fast Cals have already been done that day.

System Shutdown:

1. Remove all accessories from the cradle and store them in a cabinet during system transport.
2. At the Operator's Console, from the desktop selection area.


Desktop Selection Screen

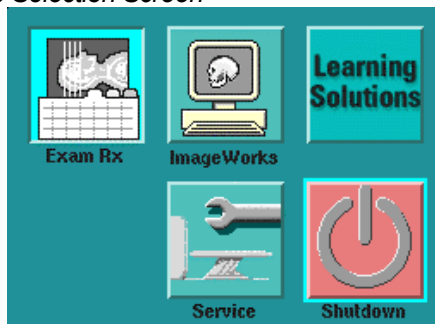

Select [Shutdown].

3. A pop-up attention box appears to begin shut down process.


Shutdown Screen

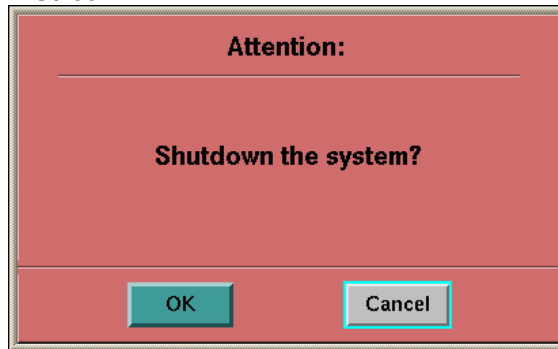

Select [OK] to begin shutdown process.
Select [Cancel] to return without starting shutdown process.

4. If [OK] is selected, a pop-up attention box appears for selection of shutdown to prepare for mobile transport or shutdown to restart the system.

Mobile Shutdown Screen

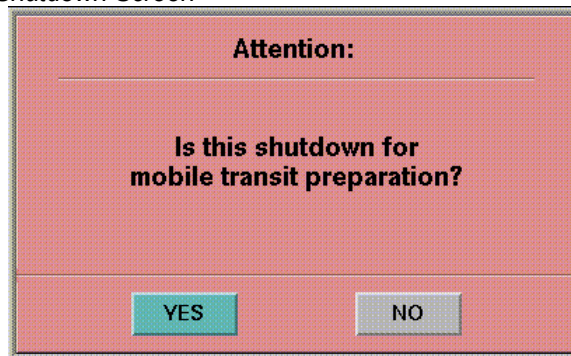

Select [Yes] to shut the system down for mobile transport.
Select [No] to shut the system down for a system restart or to power the console off.

5. A pop-up attention box is displayed instructing the user to move the gantry and table to the transport position.

Gantry Transport Setup Screen



Press the Gantry Home button on the gantry control panel. The gantry will automatically be placed in a full backward tilt position (S 30 degrees) and the table moved to full up position. The gantry and table must be set to these positions to ensure proper transport state.

6. When the system detects that it is in the proper transport position, a pop-up attention box will be displayed instructing the user to install the table transport support pole.


Table Transport Setup Screen



Select [OK] after the Table Transport Support Pole has been placed at the foot end of the table.

7. The system will automatically complete the shut down process. A series of messages will appear as the system shuts down. This will take a few minutes.
8. The system will display the following message when it is okay to turn power off at the operation console.

Restart Screen

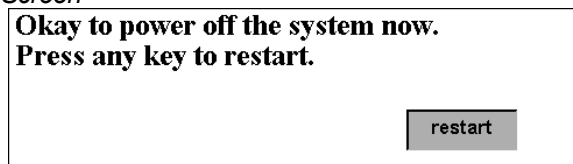


9. Verify that the red "Transport Warning Lamp" located outside the coach on the drivers side of the coach near the front of the van is not illuminated. If this light is illuminated, it means the gantry is not in the proper transport position. Restart the system and perform another shutdown for transport and ensure the gantry is tilted the full 30 degrees back.
10. Turn power off to the operator console.
11. Make sure the keyboard, mouse, trackball and chair are secured for transport with appropriate restraints.
12. Make sure the Injector is secured.
13. If your system has Cardiac Scoring option with SmartScorePro, make sure the EKG monitor is secured for transport.

Adding Names to Mobile Site List:


1. Select [Service] button in the desktop selection area.

Desktop Selection Screen

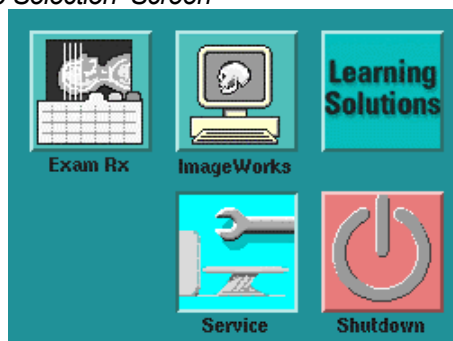


2. The Service Web Browser will appear.

Service Web Browser Screen

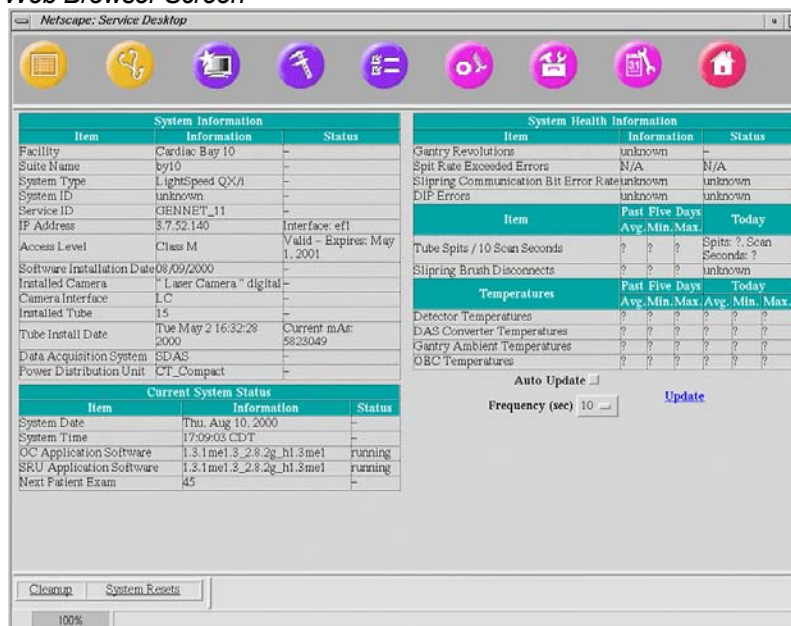

3. Select the [Utilities] button.

Service Web Utilities Screen

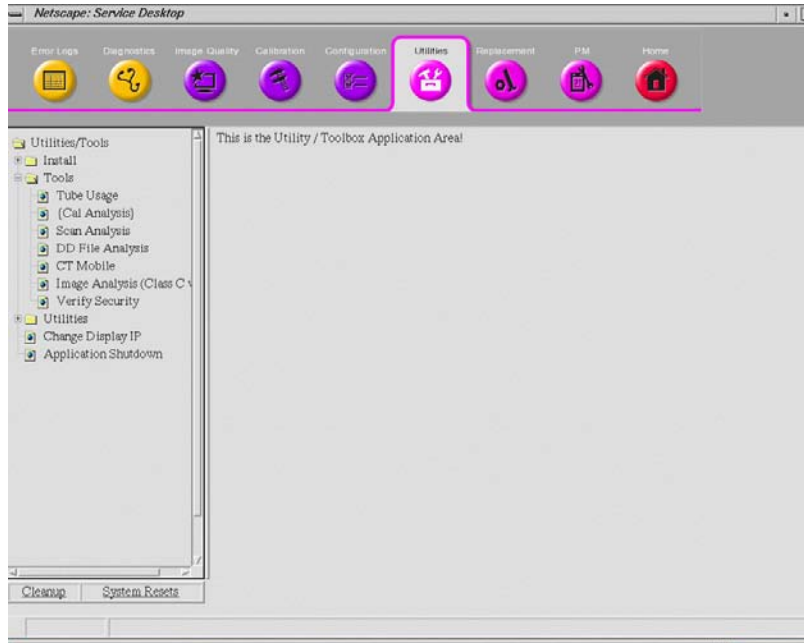

4. Select [Tools] folder.
5. Select [CT Mobile].
6. The CT Mobile Site Setup screen will be displayed.

Mobile Site Setup Screen

A screenshot of a web application titled "Mobile Site Setup" with a teal background. The main heading is "CT Mobile Site Setup". Below it, the section "Site Information:" contains several input fields: "Site Name:", "SuiteID:", "Hostname:", "IP Address:", "Netmask:", "Broadcast:", and "Default Gateway:". To the right of these fields is a large text area containing a list of hospital names: "Hospital A", "Hospital B", "Hospital C", and "Hospital F". At the bottom of the form, there are buttons: "Add", "Change", "Remove", "Search", "Select", and "Done". Below the buttons is a long, empty text input field.

- Select [Add] to add a new site to the mobile site list.

Mobile Site Setup Add Screen

Mobile Site Setup

CT Mobile Site Setup

Site Information:

Site Name:

SuiteID:

Hostname:

IP Address:

Netmask:

Broadcast:

Default Gateway:

Hospital A
Hospital B
Hospital C
Hospital F

Add Change Remove Search Select Done

ADD a new site ...

Click on [Add] and complete Site Information template.
Select [Add] to add the site.
Select [Done] to exit.

- Select [Change] to change information about a site from the mobile site list.

Mobile Site Setup Change Screen

Mobile Site Setup

CT Mobile Site Setup

Site Information:

Site Name:

SuiteID:

Hostname:

IP Address:

Netmask:

Broadcast:

Default Gateway:

Hospital A
Hospital B
Hospital C
Hospital F

Add Change Remove Search Select Done

Changing Site information ...

Select the site from list and click on [Change]. Enter corrected data in desired field for Site Information.
Click on [Change] to enter corrected information.
Select [Done] to exit.

- Select [Remove] to remove a site permanently from the mobile site list.

Mobile Site Setup Remove Screen

Select the site from the list and click on [Remove].
Select [Yes] to confirm removal of the site from the list.
Select [No] to exit removal mode without removing the site.
Select [Done] to exit.

- Select [Search] to find a site on the mobile site list by hospital name or network hostname.

Mobile Site Setup Search Screen

Select [Search] to bring up the Search pop-up.
Enter hospital name or network hostname and select [Search]. Select [Done] to exit Search screen.

7. Select [Done] to exit the CT Mobile Site Setup screen.

Mobile Site Setup Exit Screen

Select [Yes] to exit and return to the Utilities tab in the Service Web Browser.
Select [No] to cancel and return to the CT Mobile Site Setup screen.



GE Medical Systems

GE Medical Systems: Telex 3797371
P.O. Box 414, Milwaukee, Wisconsin 53201 U.S.A.
(Asia, Pacific, Latin America, North America)

GE Medical Systems — Europe: Telex 698626
283, rue de la Minière, B.P. 34, 78533 Buc Cedex
France