GE Healthcare

SmartStep LightSpeed Family

Operating Documentation

2345629-100
Revision: 4

Effectivity

The information in this manual applies to the following CT and PET/CT Scanners:

LightSpeed 2.X (4-Slice)

HiSpeed QX/i (r-Slice MDAS)

LightSpeed 3.X (4/8-Slice)

LightSpeed 4.X (16-Slice)

Discovery ST (4/8/16-Slice PET)

LightSpeed 5.X (LS16 & Pro16)

Discovery VCT (PET)

LightSpeed 7.X (VCT)

LEGAL NOTES

TRADEMARKS

Adobe, the Adobe logo, Acrobat, the Acrobat logo, Exchange, and PostScript are trademarks of Adobe Systems Incorporated or its subsidiaries and may be registered in certain jurisdictions.

Microsoft is a registered trademark and Windows is a trademark of Microsoft Corporation.

All other products and their name brands are trademarks of their respective holders.

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WARNING (EN)

- This Service Manual is available in English only.
- If a customer's service provider requires a language other than English, it is the customer's responsibility to provide translation services.
- Do not attempt to service the equipment unless this service manual has been consulted and is understood.
- Failure to heed this warning may result in injury to the service provider, operator, or patient, from electric shock or from mechanical or other hazards.

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DAMAGE IN TRANSPORTATION

All packages should be closely examined at time of delivery. If damage is apparent, have notation "damage in shipment" written on all copies of the freight or express bill before delivery is accepted or "signed for" by a General Electric representative or a hospital receiving agent. Whether noted or concealed, damage MUST be reported to the carrier immediately upon discovery, or in any event, within 14 days after receipt, and the contents and containers held for inspection by the carrier. A transportation company will not pay a claim for damage if an inspection is not requested within this 14 day period.

To file a report:

- Call 1-800-548-3366 and use option 8.
- Fill out a report on <http://us44hdd21/sctq/InstallFulfill/InstalFulfillment.htm>
- Contact your local service coordinator for more information on this process.

Rev. Jan. 5, 2005

CERTIFIED ELECTRICAL CONTRACTOR STATEMENT

All electrical Installations that are preliminary to positioning of the equipment at the site prepared for the equipment shall be performed by licensed electrical contractors. In addition, electrical feeds into the Power Distribution Unit shall be performed by licensed electrical contractors. Other connections between pieces of electrical equipment, calibrations and testing shall be performed by qualified GE Medical personnel. The products involved (and the accompanying electrical installations) are highly sophisticated, and special engineering competence is required. In performing all electrical work on these products, GE will use its own specially trained field engineers. All of GE's electrical work on these products will comply with the requirements of the applicable electrical codes.

The purchaser of GE equipment shall only utilize qualified personnel (i.e., GE's field engineers, personnel of third-party service companies with equivalent training, or licensed electricians) to perform electrical servicing on the equipment.

IMPORTANT...X-RAY PROTECTION

X-ray equipment if not properly used may cause injury. Accordingly, the instructions herein contained should be thoroughly read and understood by everyone who will use the equipment before you attempt to place this equipment in operation. The General Electric Company, Medical Systems Group, will be glad to assist and cooperate in placing this equipment in use.

Although this apparatus incorporates a high degree of protection against x-radiation other than the useful beam, no practical design of equipment can provide complete protection. Nor can any practical design compel the operator to take adequate precautions to prevent the possibility of any persons carelessly exposing themselves or others to radiation.

It is important that anyone having anything to do with x-radiation be properly trained and fully acquainted with the recommendations of the National Council on Radiation Protection and Measurements as published in NCRP Reports available from NCRP Publications, 7910 Woodmont Avenue, Room 1016, Bethesda, Maryland 20814, and of the International Commission on Radiation Protection, and take adequate steps to protect against injury.

The equipment is sold with the understanding that the General Electric Company, Medical Systems Group, its agents, and representatives have no responsibility for injury or damage which may result from improper use of the equipment.

Various protective materials and devices are available. It is urged that such materials or devices be used.

LITHIUM BATTERY CAUTIONARY STATEMENTS

CAUTION

Risk of Explosion.

Danger of explosion if battery is incorrectly replaced.

Replace only with the same or equivalent type recommended by the manufacturer. Discard used batteries according to the manufacturer's instructions.

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OMISSIONS & ERRORS

Customers, please contact your GE Sales or Service representatives.

GE personnel, please use the Healthcare PQR Process to report all omissions, errors, and defects in this publication.

Revision History

Revision	Date	Reason for change
0	06/17/2002	New Manual for SmartStep
1	06/26/2003	Added "HiSpeed QX/i" product name.
2	06/03/2004	Added LS 5.X and DST.
3	11/23/2004	PQR 13026983: Updated Prerequisites.
4	4/10/07	PQR13092807: Chapter 3: Updated information for 7.X VCT systems. Updated front and back covers, warnings section with new information and additional languages, and company name.

Table of Contents

Preface

Publication Conventions 17

Section 1.0

Safety & Hazard Information 17

1.1 Text and Character Representation..... 17

1.2 Graphical Representation 18

Section 2.0

Buttons, Switches and Keyboard Inputs (Hard & Soft Keys)..... 18

Chapter 1

SmartStep Theory..... 19

Section 1.0

Operational Overview 19

Section 2.0

Hand-Held Controller 20

Section 3.0

Prerequisites..... 22

Chapter 2

Installation - LS2.X, 3.X, 4.X & PET DST 23

Section 1.0

Overview 23

Section 2.0

Procedure..... 25

2.1 Preparation - Cover Removal 25

2.2 Installation - Faceplates..... 26

2.3 Installation - STC Connections (CT Scanners)..... 28

2.4 Installation - STC Connections (DST Scanners) 29

Chapter 3

Installation - LightSpeed 5.X and 7.X (CT and PET) 31

Section 1.0

Overview 31

Section 2.0

Procedure..... 33

2.1 Preparation - Cover Removal 33

2.2 Installation - Faceplates..... 35

2.3 Installation - Cables (LS 16, Pro16, and VCT)..... 37

Chapter 4

Finalization (All Systems)..... 39

Section 1.0

Installation - SmartStep Cradle Handle and HHC Holder 39

Section 2.0

Installation - Console Modifications 40

Section 3.0

Finishing Up - Testing & Cover Re-Installation 40

Preface

Publication Conventions

Section 1.0

Safety & Hazard Information

1.1 Text and Character Representation

Within this publication, different paragraph and character styles have been used to indicate potential hazards. Paragraph prefixes, such as hazard, caution, danger and warning, are used to identify important safety information. Text (Hazard) styles are applied to the paragraph contents that is applicable to each specific safety statement. Words describe the type of potential hazard that may be encountered and are placed immediately before the paragraph it modifies. Safety information will normally include:

- Type of potential Hazard
- Nature of potential injury
- Causative condition
- How to avoid or correct the causative condition

EXAMPLES OF HAZARD STATEMENTS USED

It is important that you read and understand hazard statements, and not just ignore them.

DANGER SEVERE PERSONAL INJURY

DANGER IS USED WHEN A HAZARD EXISTS WHICH WILL CAUSE SEVERE PERSONAL INJURY OR DEATH IF INSTRUCTIONS ARE IGNORED. THEY CAN INCLUDE:

- **ELECTROCUTION**
- **CRUSHING**
- **RADIATION**

WARNING SERIOUS PERSONAL INJURY

WARNING IS USED WHEN A HAZARD EXISTS WHICH COULD OR CAN CAUSE SERIOUS PERSONAL INJURY OR DEATH IF INSTRUCTIONS ARE IGNORED. THEY CAN INCLUDE:

- **POTENTIAL FOR SHOCK**
- **EXPOSED WIRES**
- **FAILURE TO TAG AND LOCKOUT SYSTEM POWER COULD ALLOW FOR UN-COMMAND MOTION.**

CAUTION Minor Personal Injury

Caution is used when a hazard exists which can or could cause minor injury to self or others if instructions are ignored. They include for example:

- **Loss of critical patient data**
- **Crush or pinch points**
- **Sharp objects**


NOTICE
Property
Damage
Only

Notice is used when a hazard is present that can cause property damage but has absolutely no personal injury risk. They can include:

- Disk drive will be destroyed
- Internal mechanical damage, such as to the x-ray tube
- Coasting the rotor through resonance.

1.2 Graphical Representation

Important information will always be preceded by the exclamation point  contained within a triangle, as seen throughout the Preface chapter. In addition to text, several different graphical icons (symbols) may be used to make you aware of specific types of hazards that could possibly cause harm.

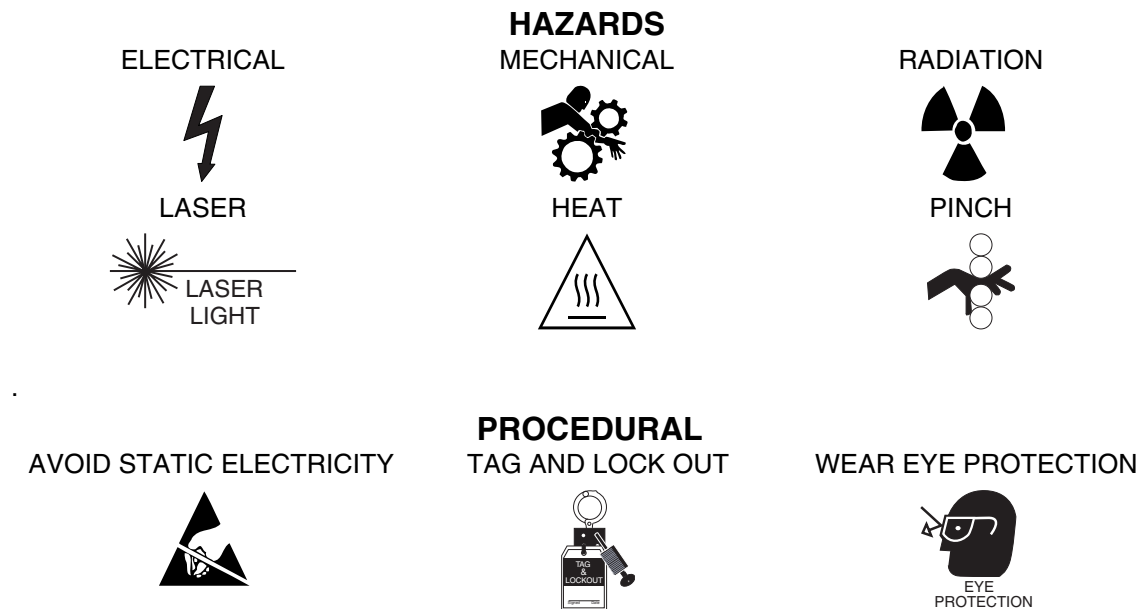Section 2.0

Buttons, Switches and Keyboard Inputs (Hard & Soft Keys)

Different character styles are used to indicate actions requiring the reader to press either a hard or soft button, switch or key. Physical hardware, such as buttons and switches, are called hard keys because they are hard wired or mechanical in nature. A keyboard or on/off switch would be a hard key. Software or computer generated buttons are called soft keys because they are software generated. Software driven menu buttons are an example of such keys. Soft and hard keys are represented differently in this publication.

Example: Hard Keys A power switch **ON/OFF** or a keyboard key like **ENTER** is indicated by applying a character style that uses both over and under-lined bold text that is bold. This is a hard key.

Example: Soft Keys Whereas the computer MENU button that you would click with your mouse or touch with your hand uses over and under-lined regular text. This is a soft key.

END OF SECTION

Chapter 1

SmartStep Theory

Section 1.0

Operational Overview

SmartStep is a mode of scanning designed to be used by the Radiologist or Physician during interventional procedures. SmartStep is accomplished by using the integrated Hand-Held Controller (HHC) and foot switch.

The available slice thicknesses of SmartStep are: 2.5 mm, 5 mm, and 7.5 mm. These are created using a 4X detector configuration. Detector rows one and two create image one. Detector rows two and three create image two. Detector rows three and four create image three. The configuration for 2.5 mm slice thickness is the 4 x 1.25 detector configuration. The configuration for 5 mm slice thickness is the 4 x 2.5 detector configuration. The configuration for 7.5 mm slice thickness is the 4 x 3.75 detector configuration.

Reconstructions for different algorithms are allowed for the same slice thickness as acquired. Only the 0.8 second and 1 second full rotation scan time retros are available for SmartStep. All algorithms are available for reconstructions.

Note: **Other desktops are unavailable during SmartStep.**

Section 2.0

Hand-Held Controller

During needle biopsies, the Radiologist or Physician controls radiation with the foot switch, and uses the hand-held control (HHC) to view images on an in-room monitor. The HHC also provides alignment light, move-to-scan and cradle micro positioning controls. The Radiologist or Physician also has the option to release HHC control of the cradle, and use the cradle handle to manually position the patient in the scan field of view.

Refer to [Figure 1-1](#) and [Table 1-1](#), on page 21 for a description of each of the HHC push button functions.

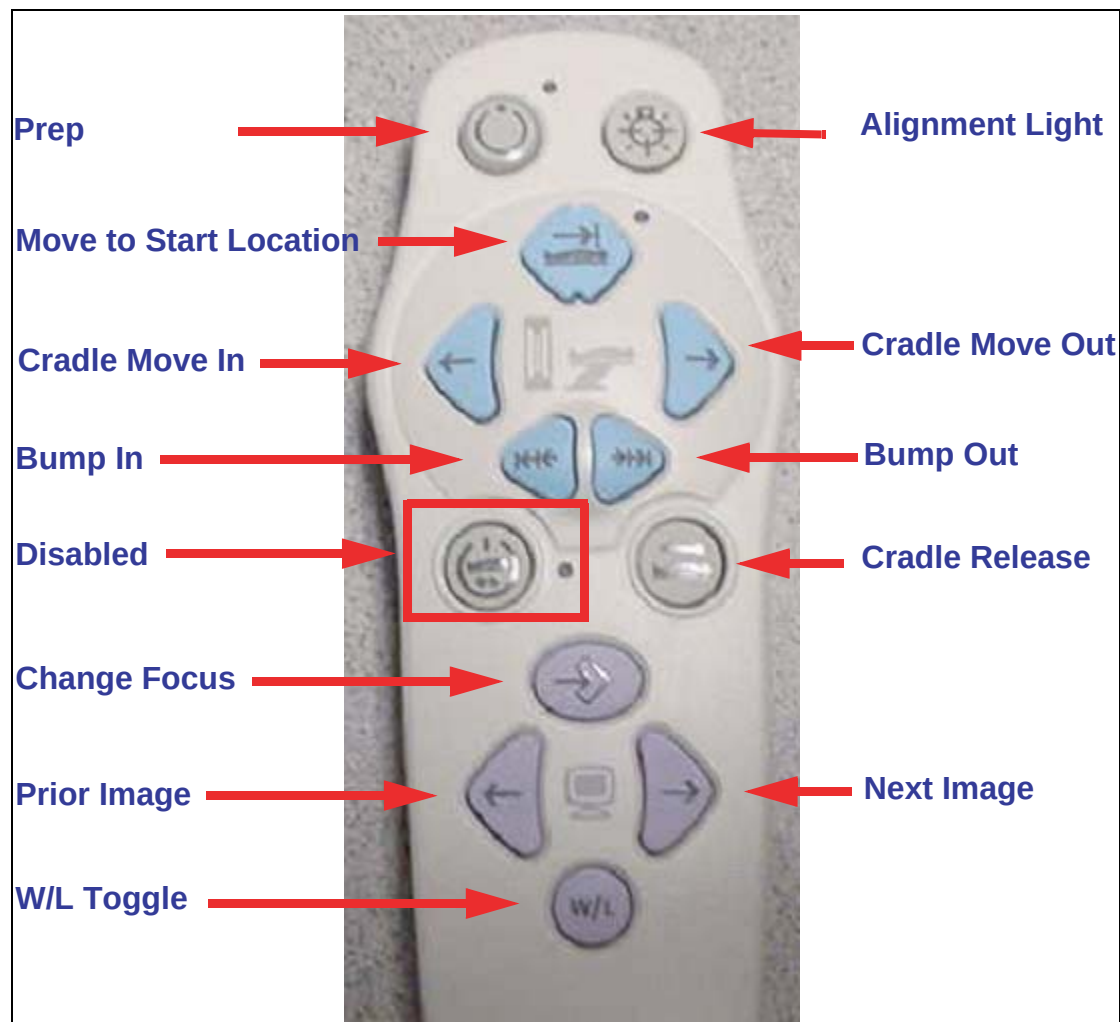
Figure 1-1 HHC Push buttons

Table 1-1 HHC Push button Functions

Button Name	Button	Description
Prep		Prepares the system for x-ray acquisitions.
Alignment Light		This button Enables the laser positioning lights.
Move to Start Location		This button positions the cradle to the start location prescribed by you in the SmartStep View/Edit screen.
Cradle Move In		The Cradle Move In button moves the cradle in towards the gantry.
Cradle Move Out		The Cradle Move In button moves the cradle away from the gantry.
Bump In		This moves the cradle a pre-defined bump distance towards the gantry. The bump distance is pre-defined on the View/Edit screen.
Bump Out		This moves the cradle a pre-defined bump distance away from the gantry. The bump distance is pre-defined on the View/Edit screen.
Disabled		This button is disabled. It is not used for SmartStep.
Cradle Release		This is a toggle button. It releases the cradle so it can be moved using the cradle handle. Pressing the button again latches the cradle.

Section 3.0

Prerequisites

- Must have all components in B7868SM Smart Step Collector - including software and hardware.
- Must have the correct hardware kit. B7710LM works for all systems other than LightSpeed 5.X, 7.X, DST HP60, and DVCT products.
- Must have B7710LN hardware kit for both PET DST HP60 and all LightSpeed 5.X and 7.X products. PET DST includes required additional PET hardware kit.

Chapter 2

Installation - LS2.X, 3.X, 4.X & PET DST

Section 1.0 Overview

The information in this Chapter applies to LightSpeed 2.X (Plus), LightSpeed 3.X (Ultra) and LightSpeed 4.X (16) CT systems, as well as Discovery ST (w/ 16/8/4-Slice CT Gantries) PET systems. The SmartStep hardware is shipped ready for immediate installation. The following installation procedure describes how to install that hardware on the gantry. See [Figure 2-1](#) for a general layout of cabling and connections. Refer to [Table 2-1, on page 24](#) for a listing of the parts in the kit.

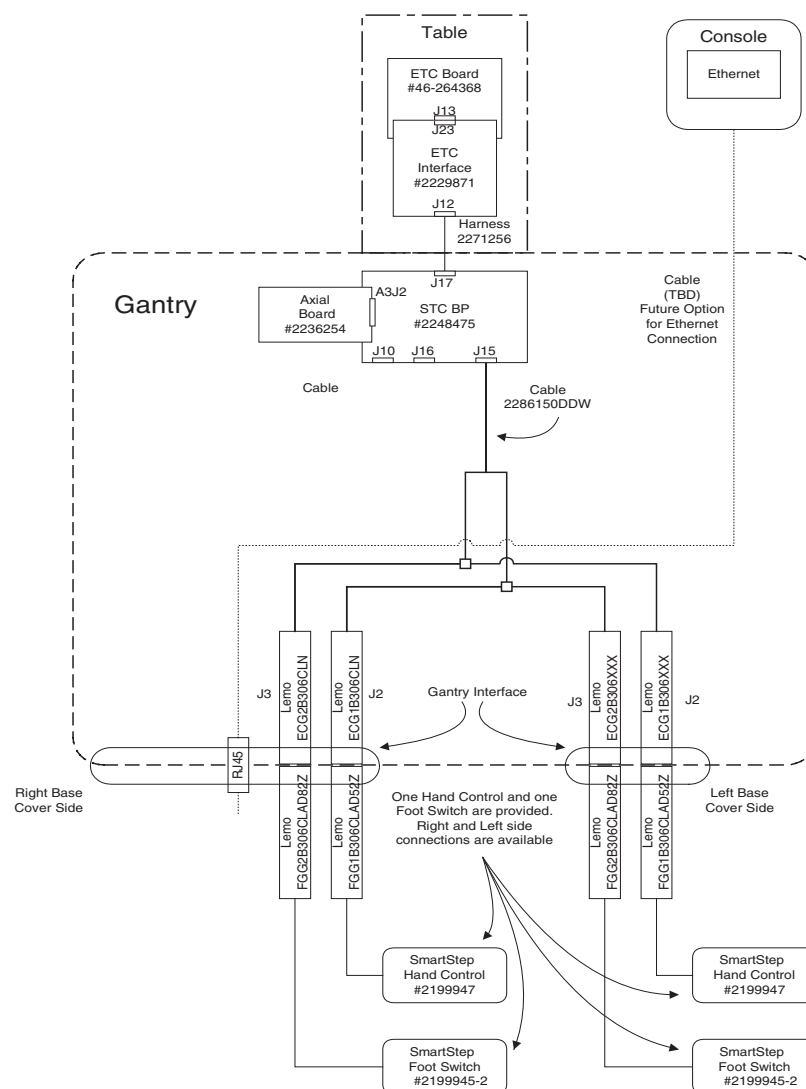
Figure 2-1 SmartStep Interconnect and Functional Drawing

See [Figure 2-2](#) for gantry's external connections.

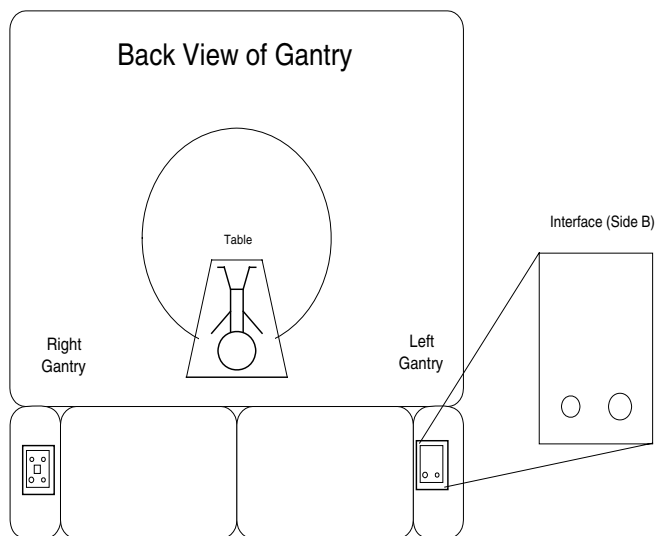
Figure 2-2 Gantry Connections

Table 2-1 2333973 SmartStep Kit Parts List (B7868WP)

Qty	Part Number	Description
1	2199947	IHHC CT Sonata hand held control (includes support bracket)
1	2199945-2	Foot switch, CT Sonata Foot Switch for H2
1	2286150	Harness, SmartStep cable - bulkhead to STC
1	2288882-2	Bracket, 8 hole bulkhead for SmartStep
1	2288882-5	Bracket, 7 hole bulkhead for SmartStep
1	2288882-3	Bracket, 11 hole bulkhead for SmartStep and Cardiac or Pulmonary
1	46-302200P10	Rating Plate
1	2206664	Cradle Handle
1	2206657	Cradle Pan and Rail Assembly
1	2206663	Hand Held Control Holder
1	2222881-2	SmartStep Software Option Key
6	46-208758P5	Tie Wraps
1	46-328417P4	Screw, Ground
1	46-328432P2	Washer, Lock, Ground
1	46-328430P2	Washer, Flat, Ground
2	2290289-2	Connector Cover, Large
2	2290289	Connector Cover, Small
2	46-328417P74	Connector Cover Screws
1	2345629-100	SmartStep - LightSpeed Plus & HiSpeed QX/i Installation Manual
1	2355108	Cable Adapter

Section 2.0 Procedure

DANGER RISK OF ELECTRIC SHOCK. USE PROPER LOCKOUT/TAGOUT PROCEDURES.

2.1 Preparation - Cover Removal

- 1.) Remove the two gantry side covers. Refer to your System Service Manual for any needed instruction.
- 2.) Turn **OFF** all 3 switches (Axial Drive, HVDC, 120 VAC) on the STC backplane.
- 3.) Remove both left and right Gantry side base covers (items 6 & 7, respectively, in [Figure 2-3](#)) and then both rear base covers (items 8 & 9, in [Figure 2-3](#)).

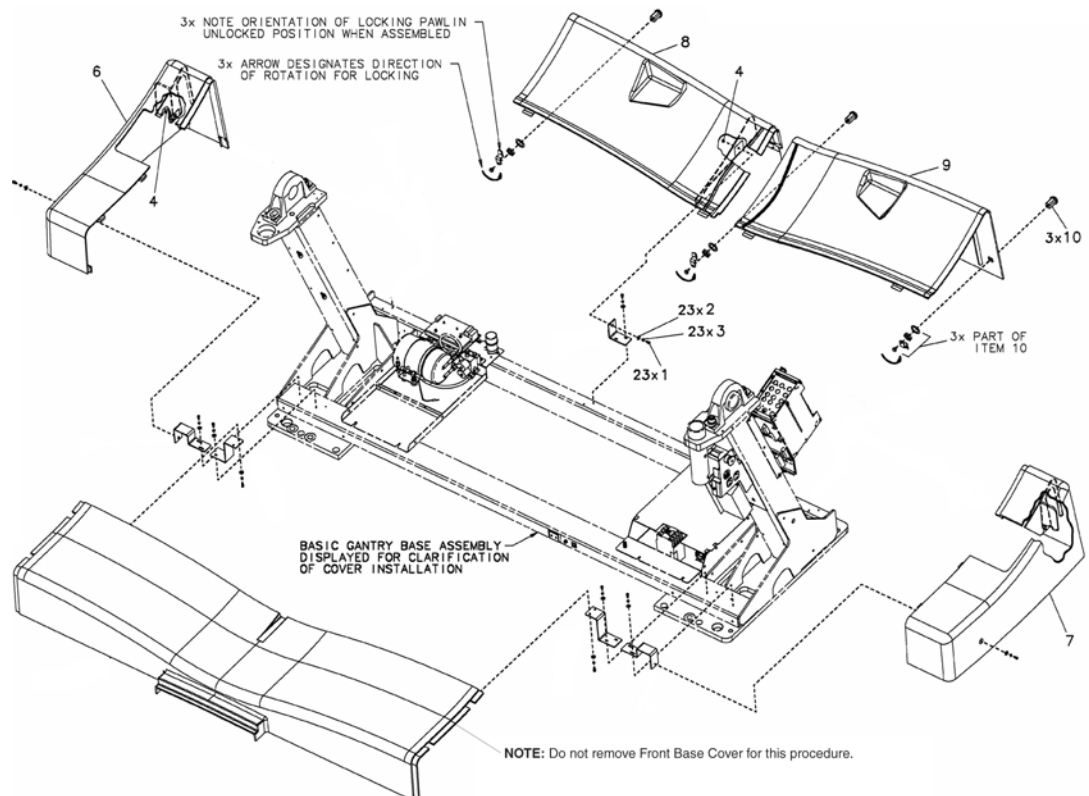
Figure 2-3 Gantry Base Covers, exploded view

Note: Do Not remove the front base cover!

2.2 Installation - Faceplates

- 1.) Remove existing faceplates - may be blanks or, if installed, Cardiac Option faceplate.
- 2.) Install the appropriate faceplate (see [Figure 2-4](#)) at indicated location per [Table 2-2](#) using the four nuts removed with the original faceplates (two per faceplate).
- 3.) Install the connector covers to the faceplates using 46-328417P74, 2290289-2, and 2290289 for each side's faceplate.

Table 2-2 Faceplate Configurations

Faceplate	Options Installed	Location
11 Hole (2288882-3)	Cardiac and SmartStep	STC Side of Gantry (Right)
8 Hole (2288882-2)	SmartStep	STC Side of Gantry (Right)
7 Hole (2288883-5)	Both Configurations	Other Side (Left)

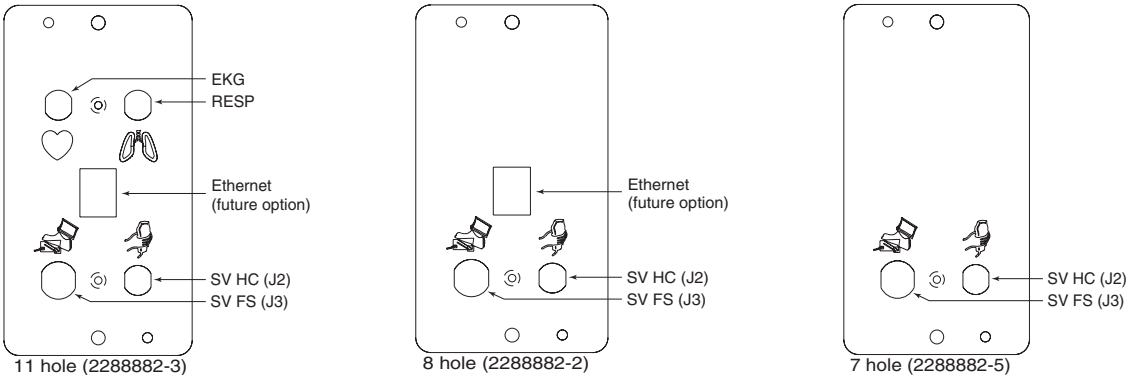
Figure 2-4 Gantry Base Faceplates

- 4.) Lay out, but do not yet mount, the included cable (2286150) along the rear cross-tube of the gantry so the connectors can reach both faceplates.
- 5.) Refer to [Figure 2-5](#). Attach the SmartStep cable to the faceplate, as follows:
 - a.) With the outer conical nut removed, insert a lemo connector through the appropriate hole in the faceplate, with the red dot facing UP (see [Figure 2-6](#), on page 27).
 - b.) Adjust the inside conical nut such that enough thread appears through the hole that the external nut can be firmly attached while keeping the connection flush with the faceplate.
 - c.) Repeat for the other three lemo connectors

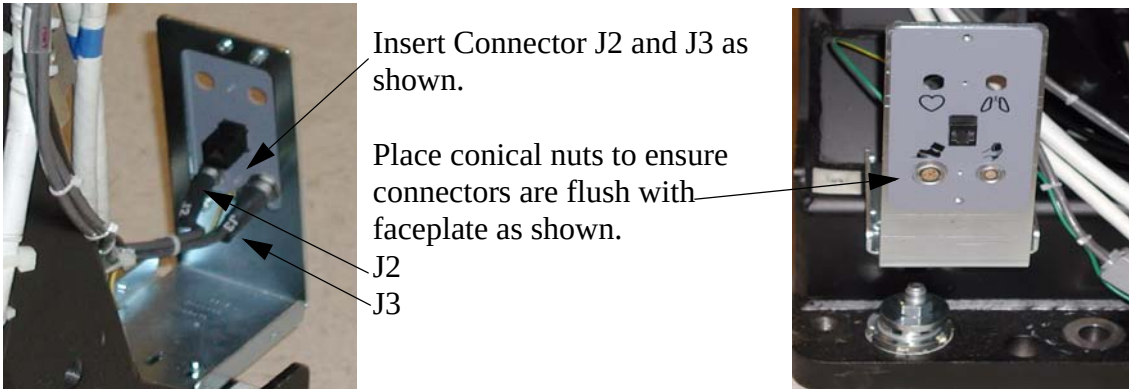
Figure 2-5 Installing the Lemo Connectors

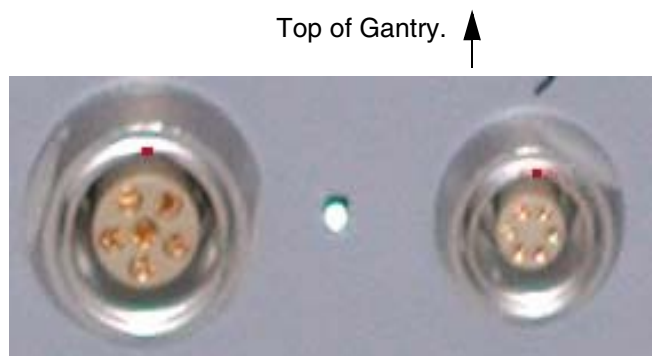
Figure 2-6 Orientation of Connectors

- 6.) Using the provided M4 screw (46-328417P4), washer (46-328430P2) and lock washer (46-328432P2), attach the SmartStep cable ground wire to the cable tie screw hole on the chassis. Remove paint as necessary to ensure a good electrical connection. See [Figure 2-7](#) for location.

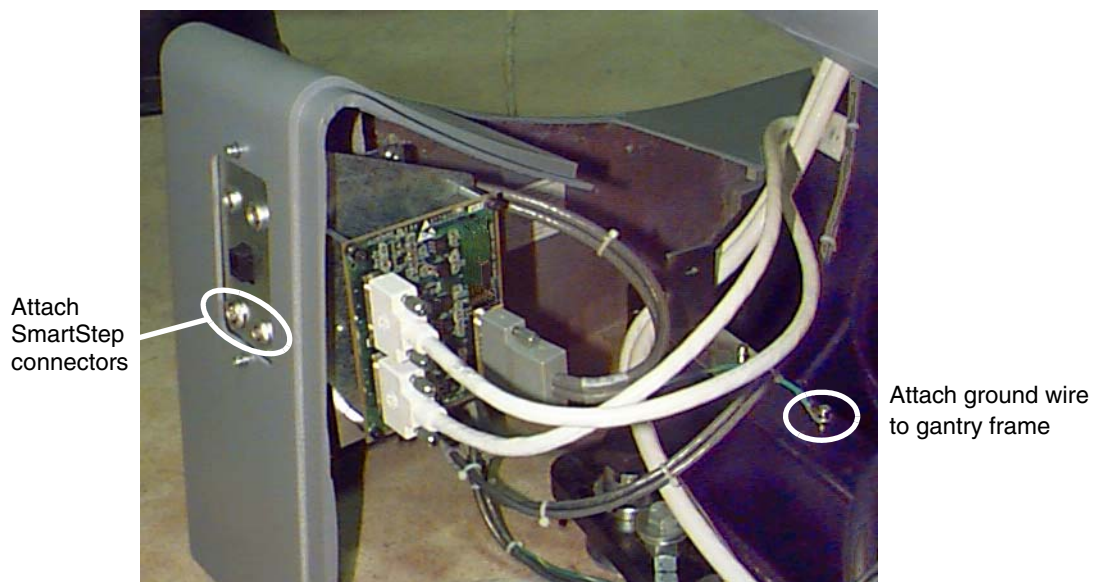
Figure 2-7 Ground Wire Installation

- 7.) Route the cables and tie wrap in place as shown in [Figure 2-8](#).

Route cable as shown along the gantry base's rear tube

Figure 2-8 SmartStep Cable Routing

2.3 Installation - STC Connections (CT Scanners)

- 1.) Connect SmartStep cable to J15 on the STC backplane. Refer to [Figure 2-9](#), for locating J15.

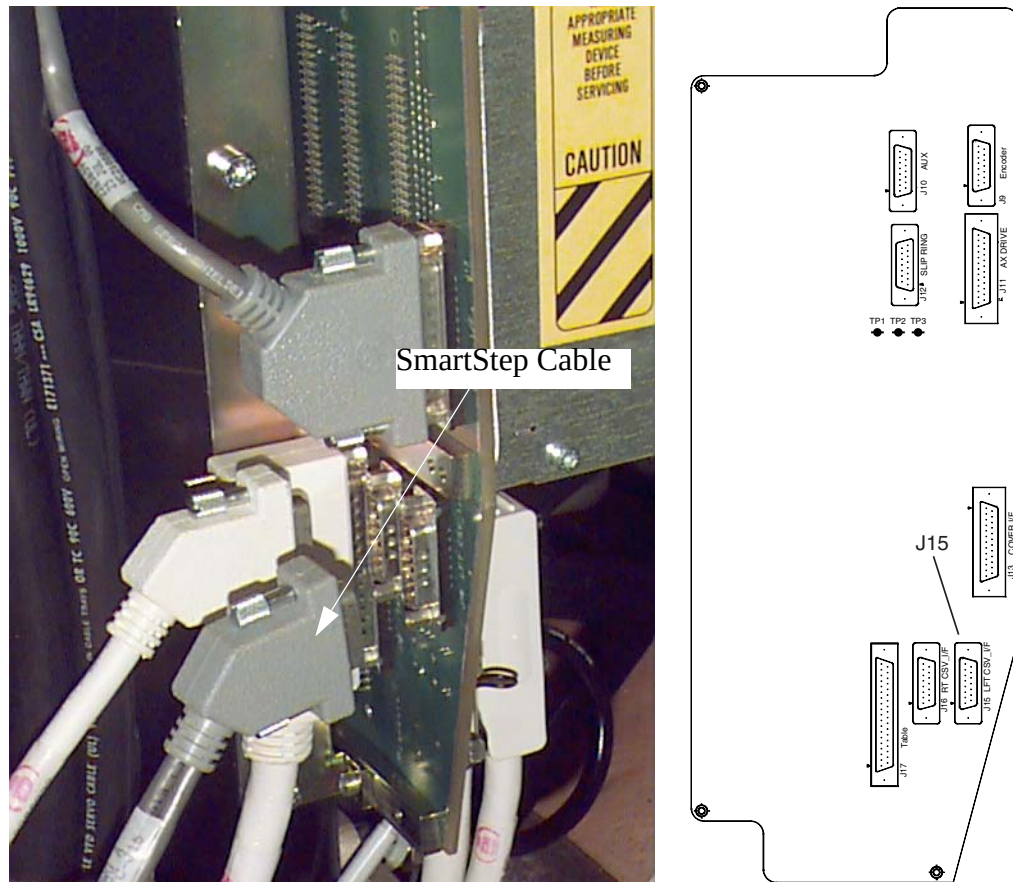
Figure 2-9 Locating connectors J15 (grey) on STC backplane

- 2.) Tie wrap the STC connection cable as shown in [Figure 2-10](#), on page 28. Be sure tie wrap is far enough down on the two cables to allow the base side cover to clear.

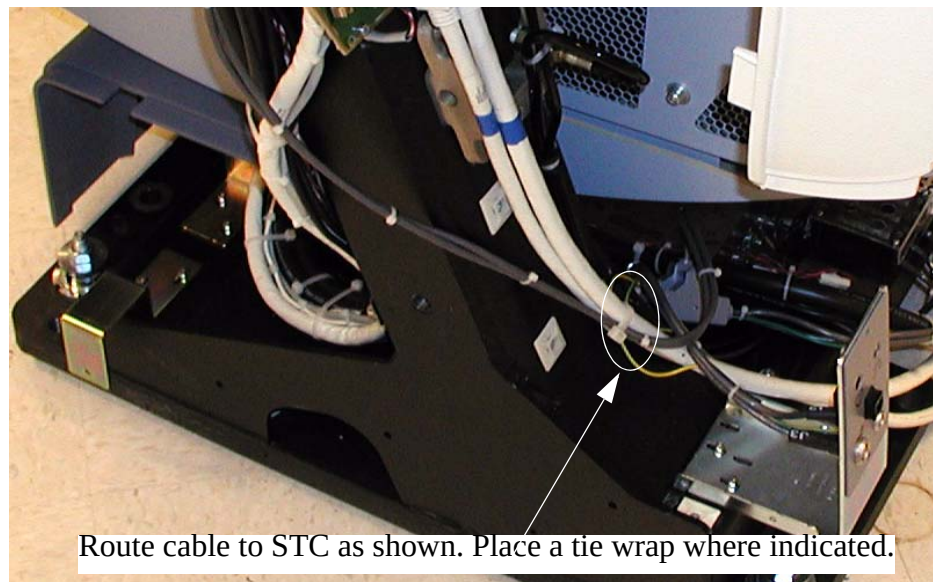
Figure 2-10 STC Cable Routing

2.4 Installation - STC Connections (DST Scanners)

Route Cable to STC as shown (Figure 2-11 - without Gating option, or Figure 2-12 - with Gating option). Place a tie-wrap, where indicated.


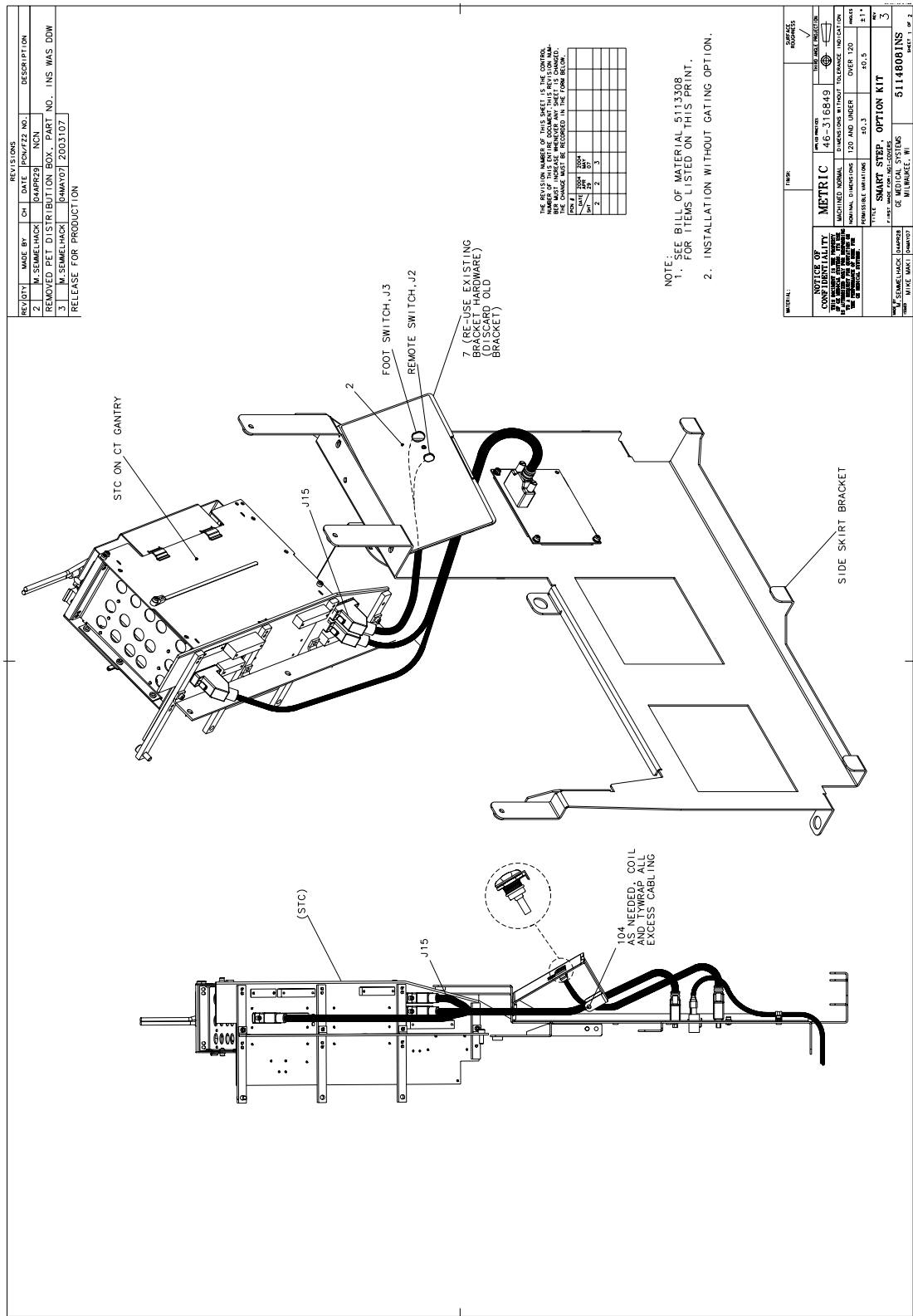
Figure 2-11 DST Cable Routing - without Gating Option

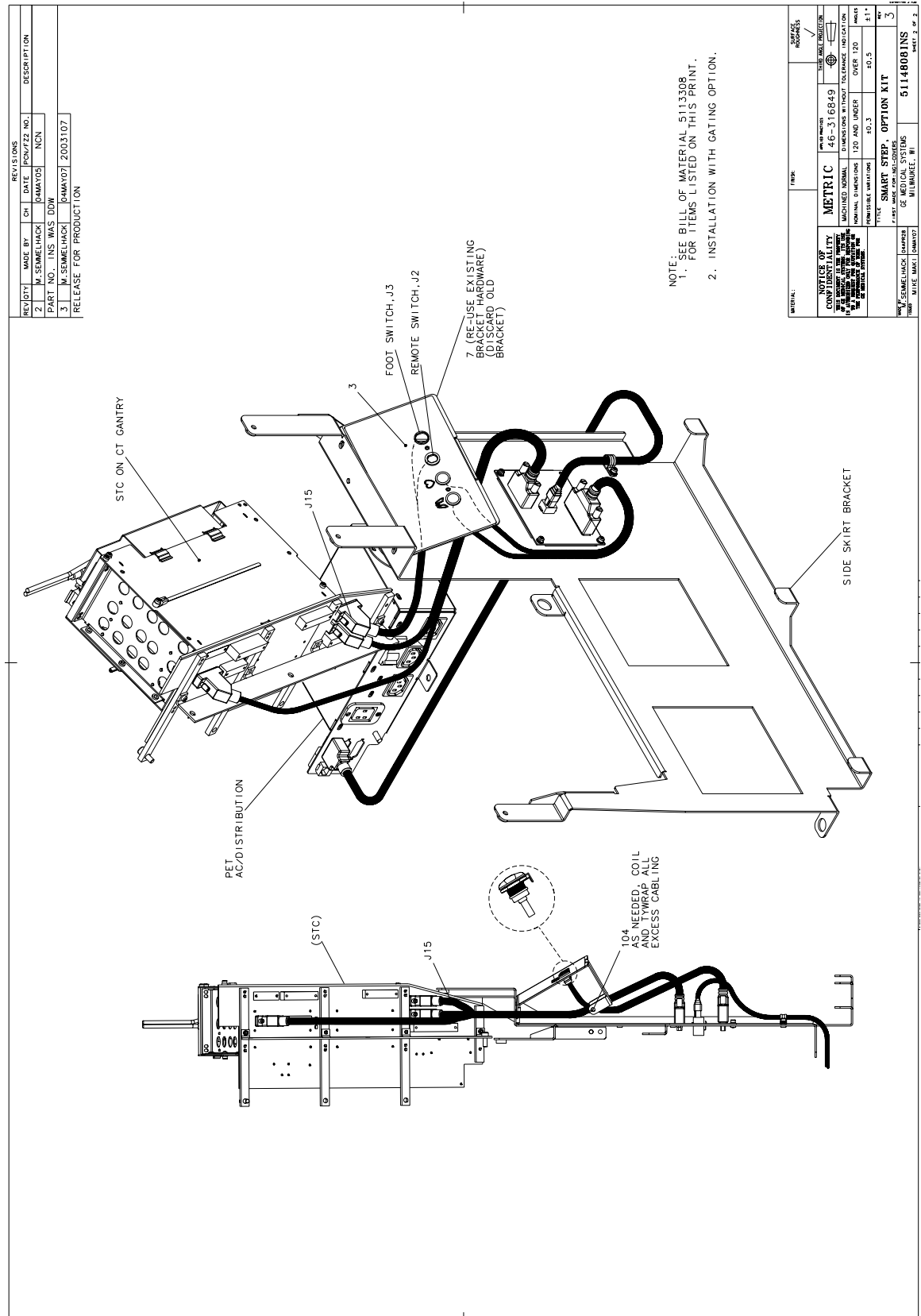Chapter 3

Installation - LightSpeed 5.X and 7.X (CT and PET)

Section 1.0 Overview

The information in this Chapter applies to LightSpeed 5.X and 7.X CT and PET systems, including LightSpeed 16, LightSpeed Pro16, and VCT systems. The information also applies to Discovery ST HP60 and DVCT systems.

The SmartStep hardware is shipped ready for immediate installation. The following installation procedure describes how to install that hardware on the gantry. See [Figure 3-1](#) for a general layout of cabling and connections. Refer to [Table 3-1](#), on page 32 for a listing of the parts in the kit.

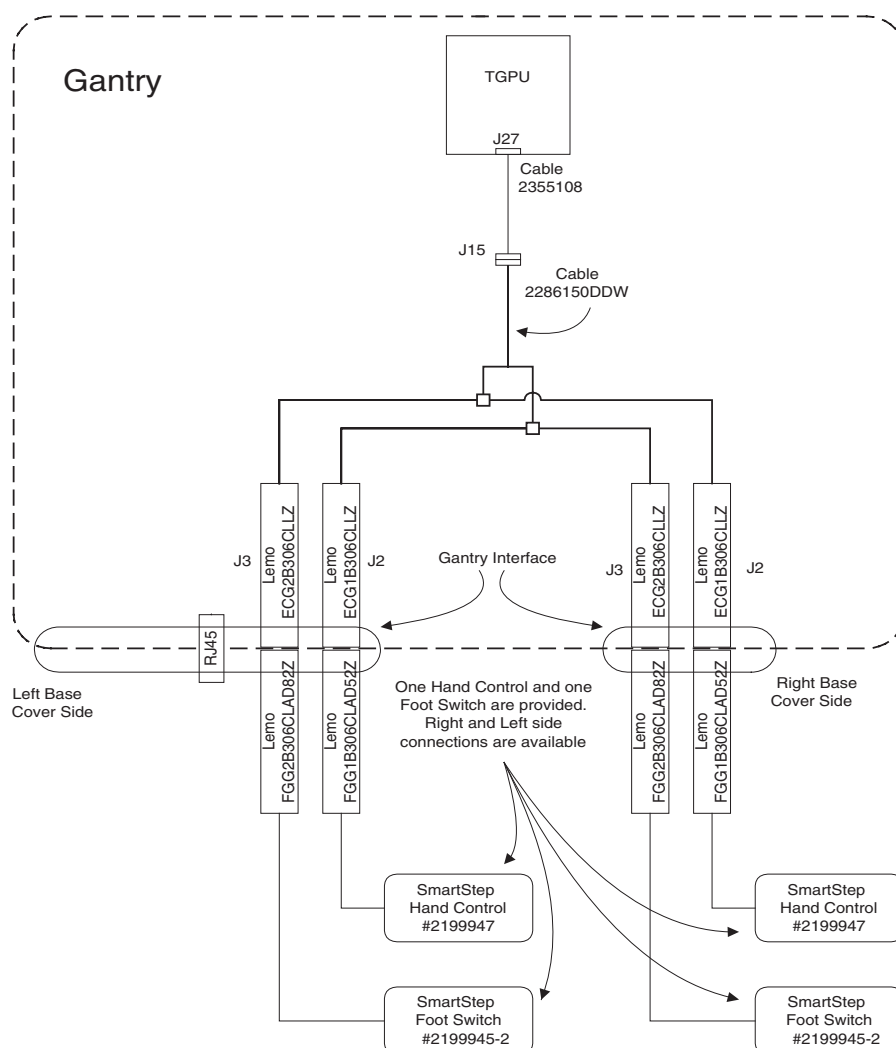
Figure 3-1 SmartStep Interconnect and Functional Drawing (LS 5.X and 7.X)

See [Figure 3-2](#) for gantry’s external connections.

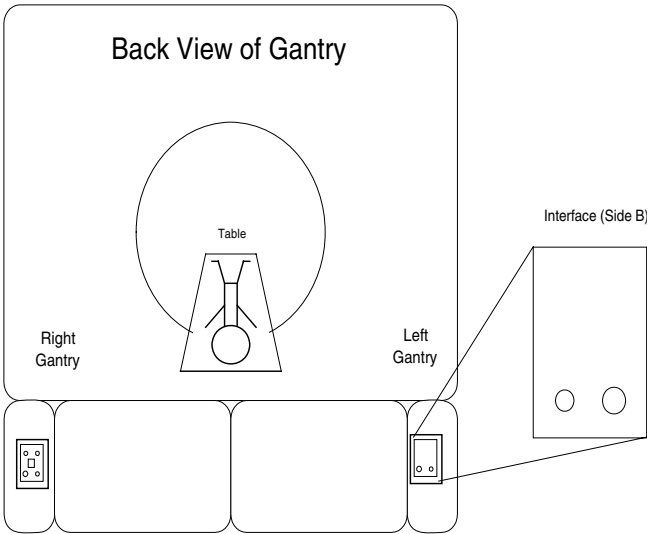
Figure 3-2 Gantry Connections

Table 3-1 2333973 SmartStep Kit Parts List (B7868SM)

Qty	Part Number	Description
1	2337414	DVD RAM Media, 9.4 GB
1	46-302233P7	Rating Plate Doc
AR	2211224-100	LightSpeed QX/I Options Installation Manual
1	2376773, incl: - 2200583-100 - 2211224-100	Boom-In-A-Room for Linux, includes: - CT/I (Sonata) Upgrade Manual - LightSpeed QX/I Options Installation Manual
1	2333973	SmartStep Kit, Controller and Footswitch

Section 2.0 Procedure

DANGER RISK OF ELECTRIC SHOCK. USE PROPER LOCKOUT/TAGOUT PROCEDURES.

2.1 Preparation - Cover Removal

- 1.) Remove the two gantry side covers. Refer to your System Service Manual for any needed instruction.
- 2.) Turn **OFF** all three gantry service switches (Axial Drive, HVDC, 120 VAC).

Figure 3-3 LightSpeed Gantry Service Switches

- 3.) Remove both left and right Gantry side base covers (items 6 & 7, respectively, in [Figure 3-4](#)) and then both rear base covers (items 8 & 9, in [Figure 3-4](#)).

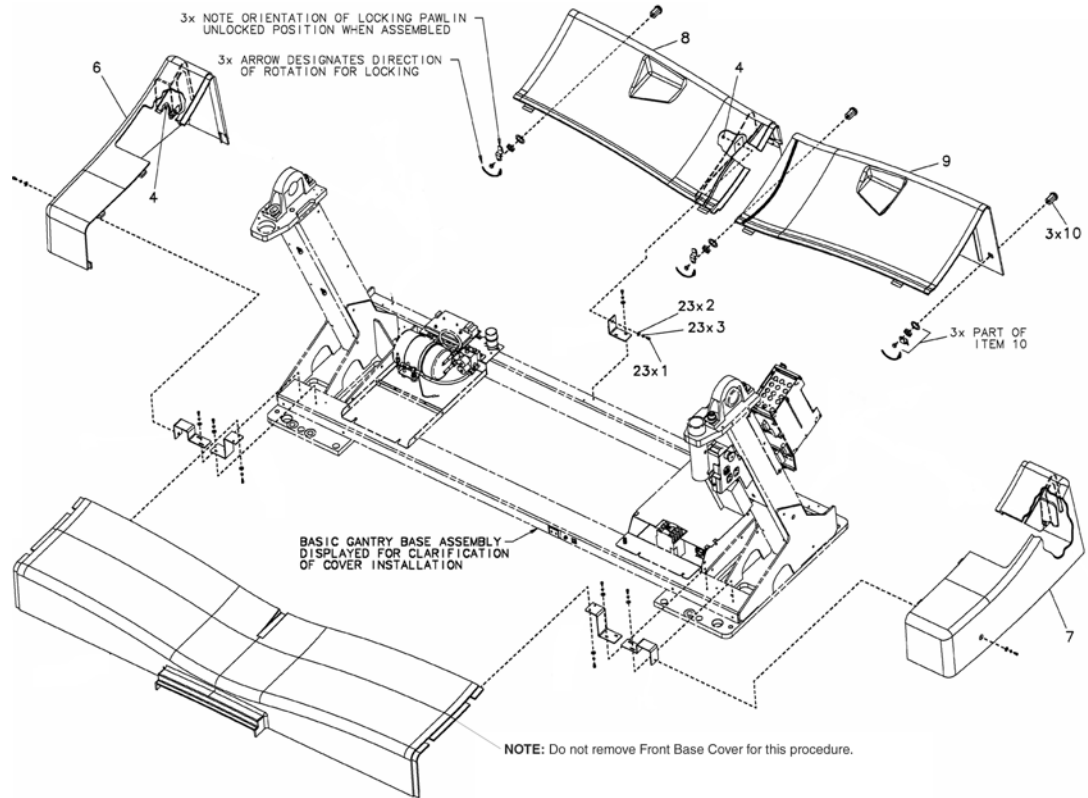
Figure 3-4 Gantry Base Covers, exploded view

Note: **Do Not remove the front base cover!**

2.2 Installation - Faceplates

- 1.) Remove existing faceplates - may be blanks or, if installed, Cardiac Option faceplate.
- 2.) Install the appropriate faceplate (see [Figure 3-5](#)) at indicated location per [Table 3-2](#) using the four nuts removed with the original faceplates (two per faceplate).

Table 3-2 Faceplate Configurations

Faceplate	Options Installed	Location
11 Hole (2288882-3)	Cardiac and SmartStep	MSUB - Left
8 Hole (2288882-2)	SmartStep	MSUB - Left
7 Hole (2288883-5)	Both Configurations	Other Side - Right

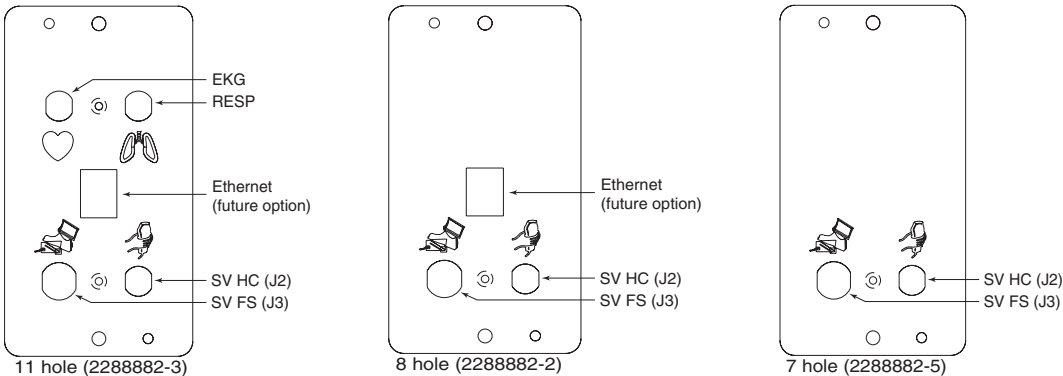
Figure 3-5 Gantry Base Faceplates

- 3.) Install the connector covers to the faceplates using 46-328417P74, 2290289-2, and 2290289 for each side's faceplate.
- 4.) Lay out, but do not yet mount, the included cable (2286150) along the rear cross-tube of the gantry so the connectors can reach both faceplates.
- 5.) Refer to [Figure 3-2, on page 32](#). Attach the SmartStep cable to the faceplate, as follows:
 - a.) With the outer conical nut removed, insert a lemo connector through the appropriate hole in the faceplate, with the red dot facing UP (see [Figure 3-6](#)).

Figure 3-6 Orientation of Connectors

- b.) Adjust the inside conical nut such that enough thread appears through the hole that the external nut can be firmly attached while keeping the connection flush with the faceplate.
- c.) Repeat for the other three lemo connectors

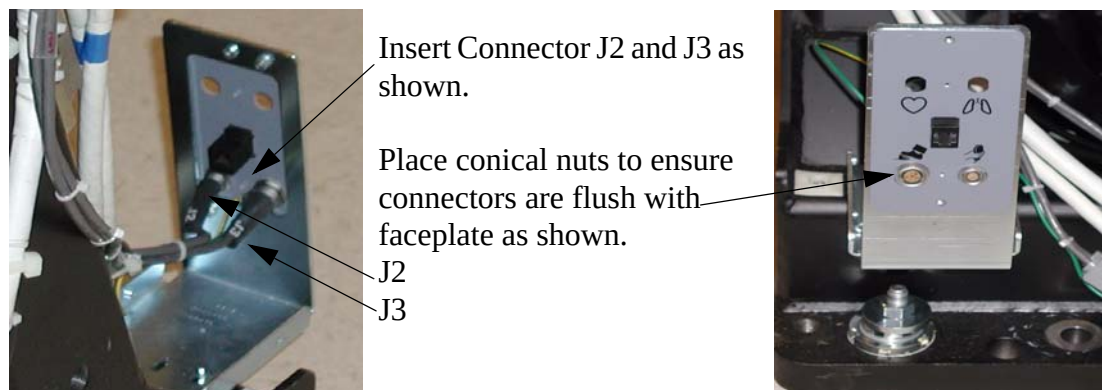
Figure 3-7 Installing the Lemo Connectors

- 6.) Using the provided M4 screw (46-328417P4), washer (46-328430P2) and lock washer (46-328432P2), attach the SmartStep cable ground wire to the cable tie screw hole on the chassis. Remove paint as necessary to ensure a good electrical connection. See [Figure 3-8](#) for location.

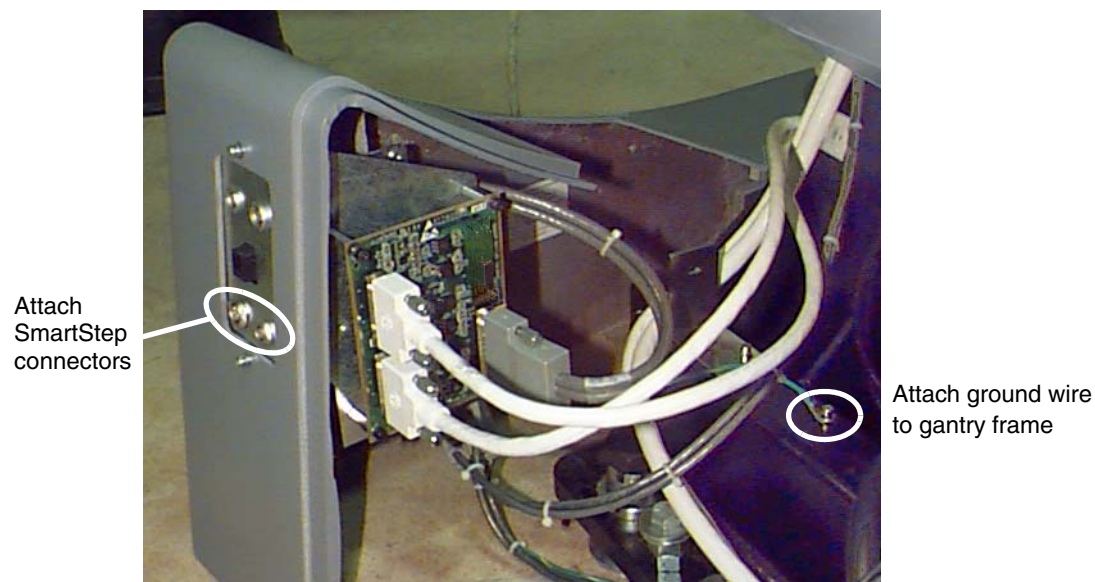
Figure 3-8 Ground Wire Installation

- 7.) Route the cables and tie wrap in place as shown in [Figure 3-9](#).

Route cable as shown along the gantry base's rear tube

Figure 3-9 SmartStep Cable Routing

2.3 Installation - Cables (LS 16, Pro16, and VCT)

- 1.) Locate the gantry MSUB. See [Figure 3-10](#).

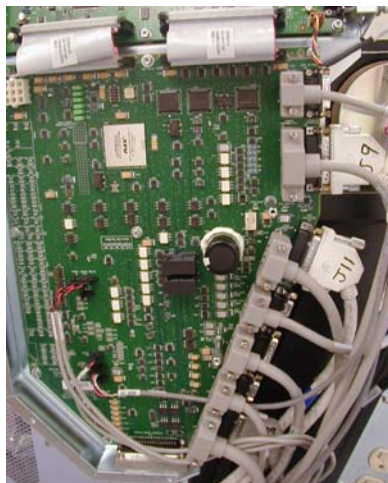
Figure 3-10 Gantry MSUB

- 2.) Locate J27. See [Figure 3-11](#).

Figure 3-11 Connector J27

- 3.) Install cable 2355108 to J27. See [Figure 3-12](#).

Figure 3-12 Install Cable 2355108

- 4.) Connect the other end to gantry cable 2286150, J15.
- 5.) Carefully route the cable along the gantry stationary frame.
- 6.) Check that covers fit around the cables.

Chapter 4

Finalization (All Systems)

Section 1.0

Installation - SmartStep Cradle Handle and HHC Holder

- 1.) Refer to [Figure 4-1](#).

Press the Cradle Release/Engage button to release HHC control.
Then, use the Cradle Handle to manually position the patient.

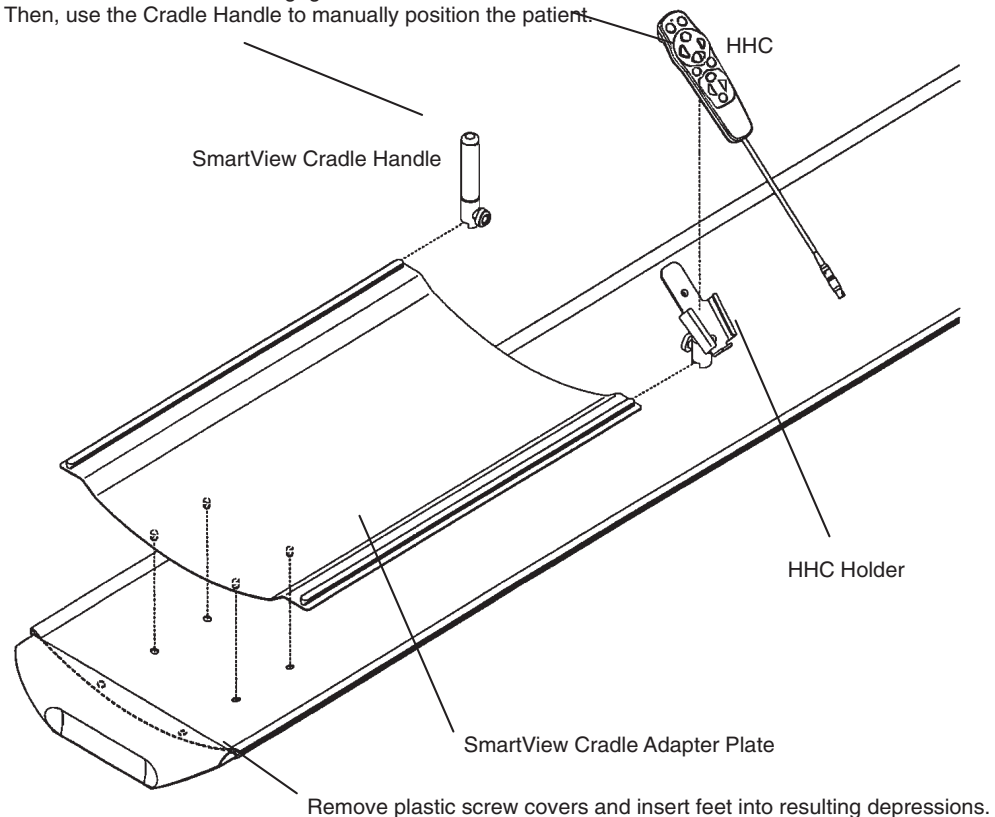
Figure 4-1 SmartStep Cradle Handle and HHC Installation

- 2.) Remove the patient pad from the cradle and set aside.
- 3.) Remove and discard the four plastic covers shipped on the cradle adaptor plate.
- 4.) Remove and discard the front four plastic caps that cover the cradle screws.
- 5.) Orient the cradle handle adaptor plate, and lay it on the cradle top so with the four “feet” on the bottom of the adaptor plate rest in the front four cradle screw holes. This prevents slippage during manual positioning.
- 6.) Return the patient pad to the cradle.
- 7.) Slide the HHC holder and cradle handle onto the cradle adaptor plate. Orient them such that the knurled knobs face away from the cradle.
- 8.) Snug the holder and handle into place by tightening the knurled knobs.

Note: The cradle handle may have some lateral slop. This is normal.

Section 2.0

Installation - Console Modifications

- 1.) Remove the modification plate from the upgrade kit.
- 2.) Refer to [Figure 4-2](#).
- 3.) Attach the new plate to the left of the existing rating plate(s).
- 4.) Take care not to cover any part of the existing plate(s). Install rating plate on console.
- 5.) Referring to your System Software Load From Cold Manual, install the SmartStep Option Key.

Figure 4-2 Operator Console Rating Plate Attachment (Left to Right: H3, GOC1 & GRE)

Section 3.0

Finishing Up - Testing & Cover Re-Installation

- 1.) Be certain all obstructions to gantry motion are clear and turn **ON** all 3 switches (Axial Drive, HVDC, 120 VAC) on the STC backplane.
- 2.) Run a system functional check, including a test scan, before installing the covers.
- 3.) Check SmartStep Operation - refer to SmartStep section of the Operator's Manual for procedure.
- 4.) Turn **OFF** all 3 switches (Axial Drive, HVDC, 120 VAC) on the STC backplane.
- 5.) Reinstall both rear base covers (See [Figure 2-3](#), on page 25 or [Figure 3-4](#), on page 34).
- 6.) Reinstall the gantry side covers. Refer to your System Service Manual for any needed instruction.

(c) 2002, General Electric Company.
GE Medical Systems, a General Electric Company, going to market as GE Healthcare.
3000 N. Grandview Boulevard
Waukesha, Wisconsin 53188
USA

www.gehealthcare.com