

MI&CT ESERVICE NOTE

SN 06121211a

Class M information is Property of GE. Disclosure to Third Parties prohibited.

Date: December 11, 2006

Subject: Lexan Removal

Effectivity: HP-60, -80, -100 systems, RT, VCT

Details: Over time the Lexan covering the collimator can get brittle and break apart leading to image quality problems. To prevent this from happening, removal of the Lexan is required.

Note: This procedure is only to be performed in conjunction with a Tube change, and only needs to be performed one time.

Resolution: See following instructions.

Personnel Requirements

Required Persons	Time
1	1 hr

Tools and Test Equipment

Item	Qty
Standard FE Toolkit	1
Loctite 242 (0.5 ml)	1

Consumables

Item	Qty
Loctite 242 (0.5 ml)	1

Safety



NOTICE


This procedure involves handling lead. Wear your lead PPE and follow lead safety handling procedure.

Prodecure

1. Lexan Removal


1. Locate the interposer plate. (See Illustration 1.)

ILLUSTRATION 1 – INTERPOSER PLATE

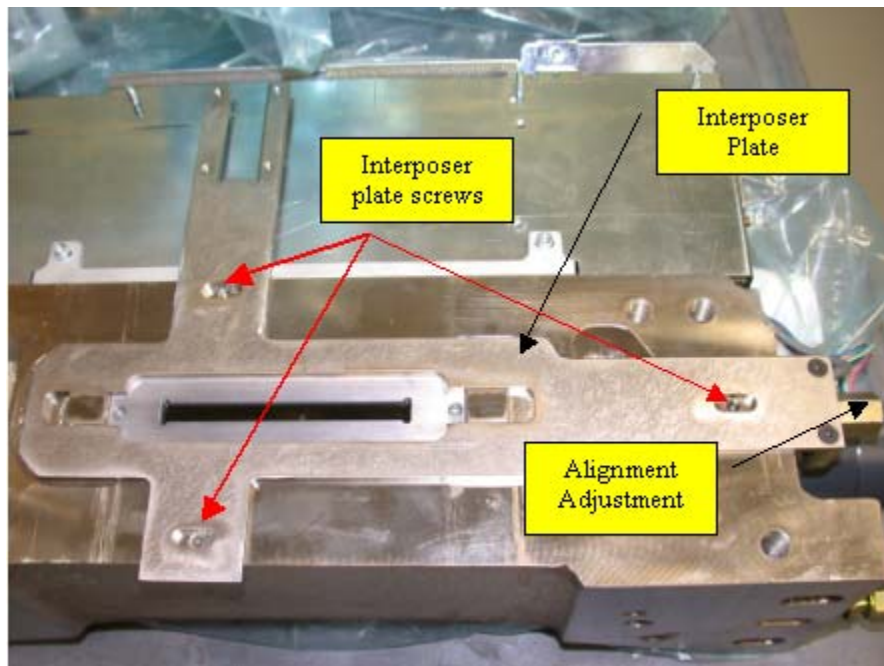

2. Before removing the interposer plate, outline its position with a marker (Illustration 2). This facilitates returning the collimator to its original position later.


ILLUSTRATION 2 – OUTLINING THE INTERPOSER PLATE

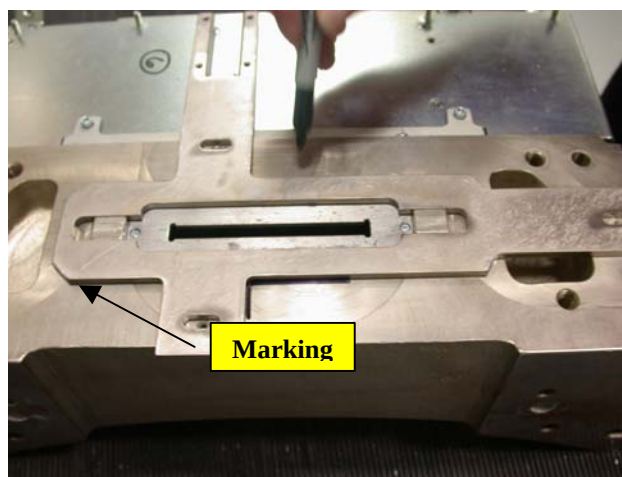

GE Healthcare MI&CT Class M eService Notes are protected by copyright (c) to General Electric Company Inc with all rights reserved. Class M eService Notes issued by GE Healthcare are proprietary to GE Healthcare Systems. Disclosure to persons not employed by GE Healthcare prohibited.
Duplication and distribution of any class of MI&CT eService Notes, in any form by entities other than GE Healthcare is strictly prohibited.


3. On VCT systems, remove 3 screws from the interposer plate as shown in Illustration 3. Set aside bracket the three screws were removed from. Do not adjust the large alignment screw, as this will help maintain interposer plate position.

ILLUSTRATION 3 – VCT INTERPOSER PLATE ALIGNMENT SCREWS

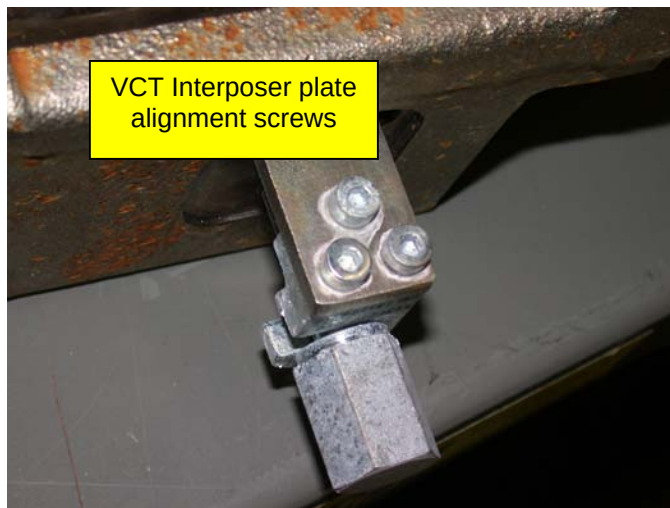

4. Drive the collimator filters out of the Field of View (FOV) to the home position (toward the Collimator Control Board) as shown in Illustration 4. This reduces the chance of contamination and potential image quality problems.


ILLUSTRATION 4 – DRIVING COLLIMATOR FILTERS TO HOME POSITION

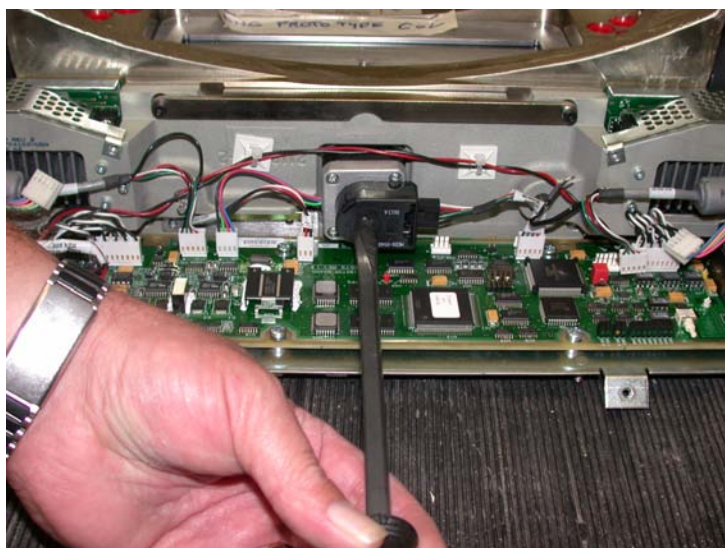

5. Remove screws from the lead aperture as shown in Illustration 5.

NOTE: VCT systems have 6 screws.

ILLUSTRATION 5 – LEAD APERTURE SCREWS

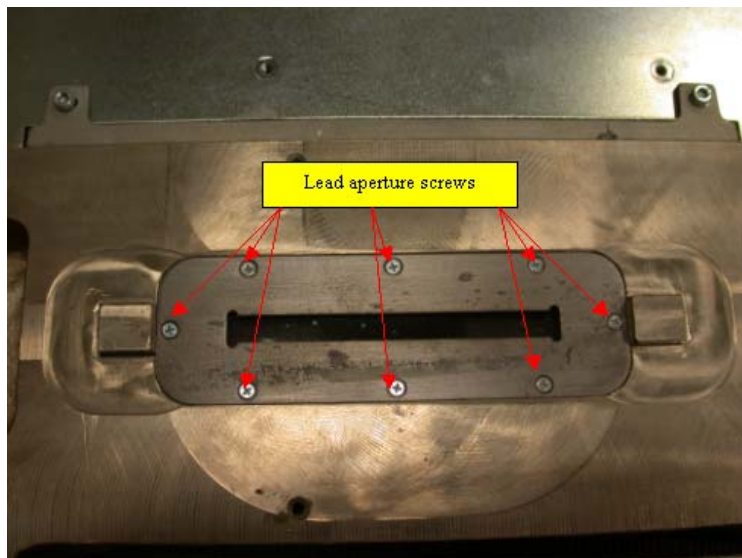

6. With a flat blade screwdriver, carefully pry out the lead aperture.
7. Clean any debris away from filter and lead aperture.
8. Remove Lexan from lead aperture and clean the underside of the aperture of any debris. (See Illustration 6.)

ILLUSTRATION 6 – CLEANING UNDERSIDE OF APERTURE



GE Healthcare MI&CT Class M eService Notes are protected by copyright (c) to General Electric Company Inc with all rights reserved. Class M eService Notes issued by GE Healthcare are proprietary to GE Healthcare Systems. Disclosure to persons not employed by GE Healthcare prohibited.
Duplication and distribution of any class of MI&CT eService Notes, in any form by entities other than GE Healthcare is strictly prohibited.

2. Reassemble Collimator

1. Using a single drop of Loctite 242 on each aperture screw, re-attach the lead aperture.
2. Re-attach the interposer plate, re-aligning it to the previously marked position. (See Illustration 2.)

Michael L. Bauer
Development Engineer – MI&CT IB Quality
MI&CT Engineering Manager Signature_____