

GE Healthcare MI&CT Class C eService Notes are copyright (c) protected to General Electric Company Inc with all rights reserved. Class C eService Notes issued by GE Healthcare are proprietary to GE Healthcare Systems. Disclosure to persons not employed by GE Healthcare or not in possession of a Advanced Service Agreement with GE Healthcare is prohibited.

Duplication and distribution of any class of eService Notes, in any form by entities other than GE Healthcare is strictly prohibited.

MI&CT ESERVICE NOTE SN07110612b

Class C information is Property of GE.

This short-term informational note is invalid 6 months after the issue date.

Date: **Reissued July 8, 2008** Original release date November 11, 2007

Subject: Hard Disk Drive Compatibility and Jumper configuration

Effectivity: ALL H1, H2, H3, GOC 1, GOC 2, and Lightspeed QX/I Consoles

System Alias and Scanner Type	Console P/Ns affected	Affiliated Software affected	Operating System P/N affected
Lightspeed QX/I, H1 H1.2, and H1.3	2180551	10.5.1_2.8.21 H1.3 M4 Software	2240216-2, 2240213-2, 2240215-2
Lightspeed QX/I, H1 H1.2, and H1.3	2180551-2	10.5.1_2.8.21 H1.3 M4 Software	2240216-2, 2240213-2, 2240215-2
Lightspeed QX/I, H1 H1.2, and H1.3	2180551-3	10.5.1_2.8.21 H1.3 M4 Software	2240216-2, 2240213-2, 2240215-2
IRIX LS2002 H2, H3 (Warp 3) and H4	2266832	All IRIX LS2005M5 Software	2356015, 2356016, 2356017
IRIX LS2002 H2, H3 (Warp 3) and H4	2266832-2	All IRIX LS2005M5 Software	2356015, 2356016, 2356017
IRIX LS2002 H2, H3 (Warp 3) and H4	2266832-3	All IRIX LS2005M5 Software	2356015, 2356016, 2356017
IRIX LS2002 H2, H3 (Warp 3) and H4	2266832-4	All IRIX LS2005M5 Software	2356015, 2356016, 2356017
IRIX LS2002 H2, H3 (Warp 3) and H4	2304732	All IRIX LS2005M5 Software	2356015, 2356016, 2356017
IRIX LS2002 H2, H3 (Warp 3) and H4	2304732-2	All IRIX LS2005M5 Software	2356015, 2356016, 2356017
IRIX LS2002 H2, H3 (Warp 3) and H4	2304732-3	All IRIX LS2005M5 Software	2356015, 2356016, 2356017
IRIX LS2002 H2, H3 (Warp 3) and H4	2304732-4	All IRIX LS2005M5 Software	2356015, 2356016, 2356017
IRIX LS2002 H2, H3 (Warp 3) and H4	2341104	All IRIX LS2005M5 Software	2356015, 2356016, 2356017
IRIX LS2002 H2, H3 (Warp 3) and H4	2341104-2	All IRIX LS2005M5 Software	2356015, 2356016, 2356017
IRIX LS2002 H2, H3 (Warp 3) and H4	2341923	All IRIX LS2005M5 Software	2356015, 2356016, 2356017
GOC 1 & 2	2315653	LINUX 1.7.10	2399296

Details:

This Service Note will discuss two (2) issues related to Hard Disk Drives in the H1, H2, H3 GOC1, GOC2 and Lightspeed QXi (SGI Octane based) consoles.

1. Operating System Hard Disk Drives (2169940-X, 9 GB OS Drive)
Scan Data Hard Disk Drive (2225004-X, 73GB 15000 RPM Scan Data Drive)

Due to obsolescence (End of Life) with Operating System and Scan Data hard disk drives, numerous types of hard disk drives exist in the installed base. A compatibility issue has been identified when certain combinations of these drives are installed in a system. The basic problem with compatibility is related to the partition size of the Scan Data drive verses the partition size of the Operating System disk drive. If a large (73 GB) Scan Data disk drive is combined in a system with a small (9 GB) Operating System disk drive, the system will produce errors related to the Operating System drive becoming full with scan files and system logs.

2. Scan Data Hard Disk Drive (2225004-9 - 73GB 15000 RPM Scan Data Drive)

A new Scan Data hard disk drive added to the FRU stock in recent months have different jumper configurations compared with previous hard disk drives used in these consoles. Though correctly configured from the vendor, concern has been raised that the jumper configuration is not shown in service documentation.

Resolution:

1. Hard Disk Drive Compatibility

If the console contains a 9 GB OS disk drive and a Scan Data disk drive needs to be replaced, order the following:

2225004-10 (73 GB destroyed to 18GB, 15000 RPM, Scan Data Drive)

If the console contains an 18, 36 or 73 GB OS disk drive and a Scan Data disk drive needs to be replaced, order the following:

2225004-9 (73 GB, 15000 RPM, Scan Data Drive)

Reverse Scenario

If the console contains an 18 or 36 GB Scan Data disk drive and an OS disk drive needs to be replaced, order the following:

2169940-72UU (9 GB, 10000 RPM, OS Disk Drive)

2169940-78UU (18 GB, 10000 RPM, OS Disk Drive)

If the console contains a 73 GB Scan Data disk drive and an OS disk drive needs to be replaced, order the following:

2169940-78UU (18 GB, 10000 RPM, OS Disk Drive)


2. Hard Disk Drive Jumper Configuration

The 2225004-9 and -10 Scan Data hard disk drives are set up as a Field Replaceable Units and will come shipped with jumper settings correctly installed as *SCSI ID 1*, *Force Single Ended* and *Term Power On*. The jumpers will not need to be altered unless you are trying to install the drive in a 16-slice system. This system has two Scan Data Hard Drives. Thus, one must have SCSI ID 1 while the other drive must have SCSI ID 2. It does not matter which drive is ID 1 and which drive is ID 2.

NOTE: Depending on hardware configuration, some 16 Slice systems with dual Scan Disk Drives may require the removal of the Forced Single-ended Bus Mode jumper.

For reference information, here are pictures detailing the default settings and locations of the jumpers on the 2225004-9 and -10 hard disk drives as well as the current 2225004-8 hard disk drive.


2225004-8 Hard Drive

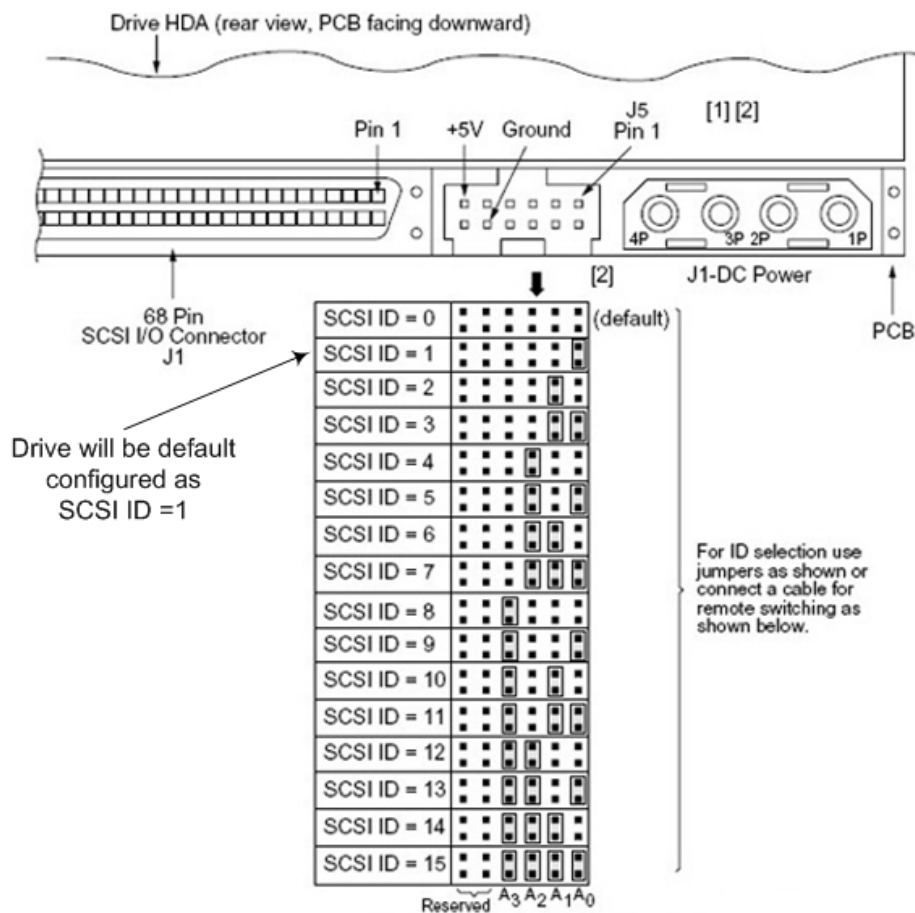

GE Healthcare MI&CT Class C eService Notes are copyright (c) protected to General Electric Company Inc with all rights reserved. Class C eService Notes issued by GE Healthcare are proprietary to GE Healthcare Systems. Disclosure to persons not employed by GE Healthcare or not in possession of a Advanced Service Agreement with GE Healthcare is prohibited.

Duplication and distribution of any class of eService Notes, in any form by entities other than GE Healthcare is strictly prohibited.


2225004-8 Hard Drive

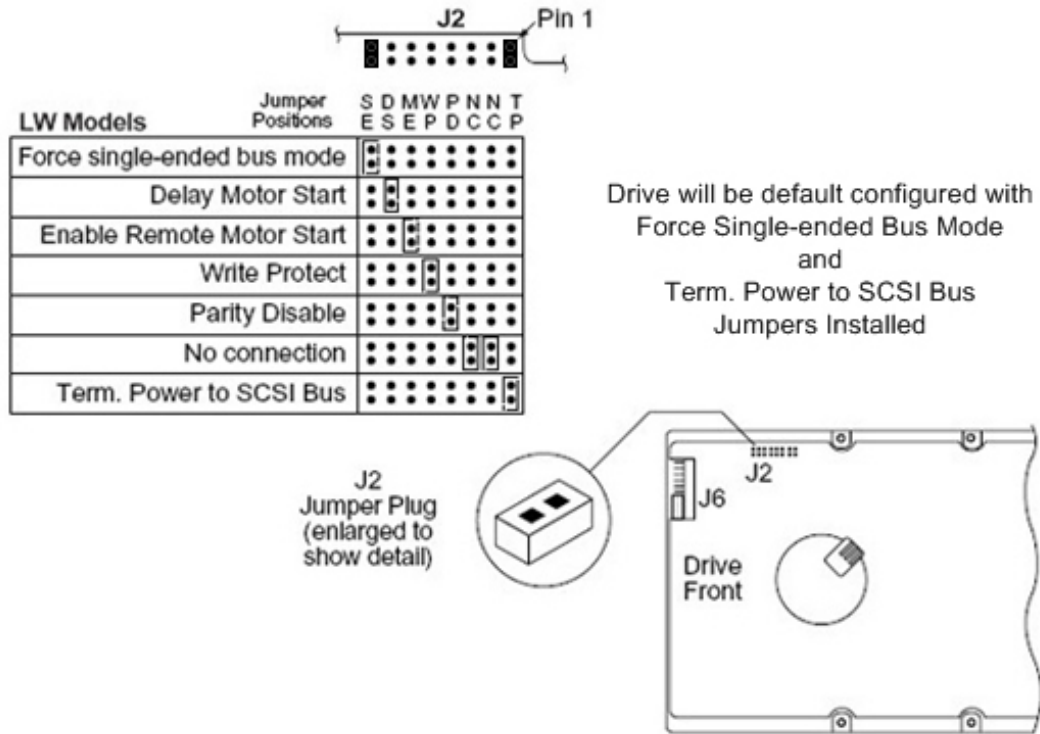

GE Healthcare MI&CT Class C eService Notes are copyright (c) protected to General Electric Company Inc with all rights reserved. Class C eService Notes issued by GE Healthcare are proprietary to GE Healthcare Systems. Disclosure to persons not employed by GE Healthcare or not in possession of a Advanced Service Agreement with GE Healthcare is prohibited.

Duplication and distribution of any class of eService Notes, in any form by entities other than GE Healthcare is strictly prohibited.


2225004-9 & -10 Hard Drives

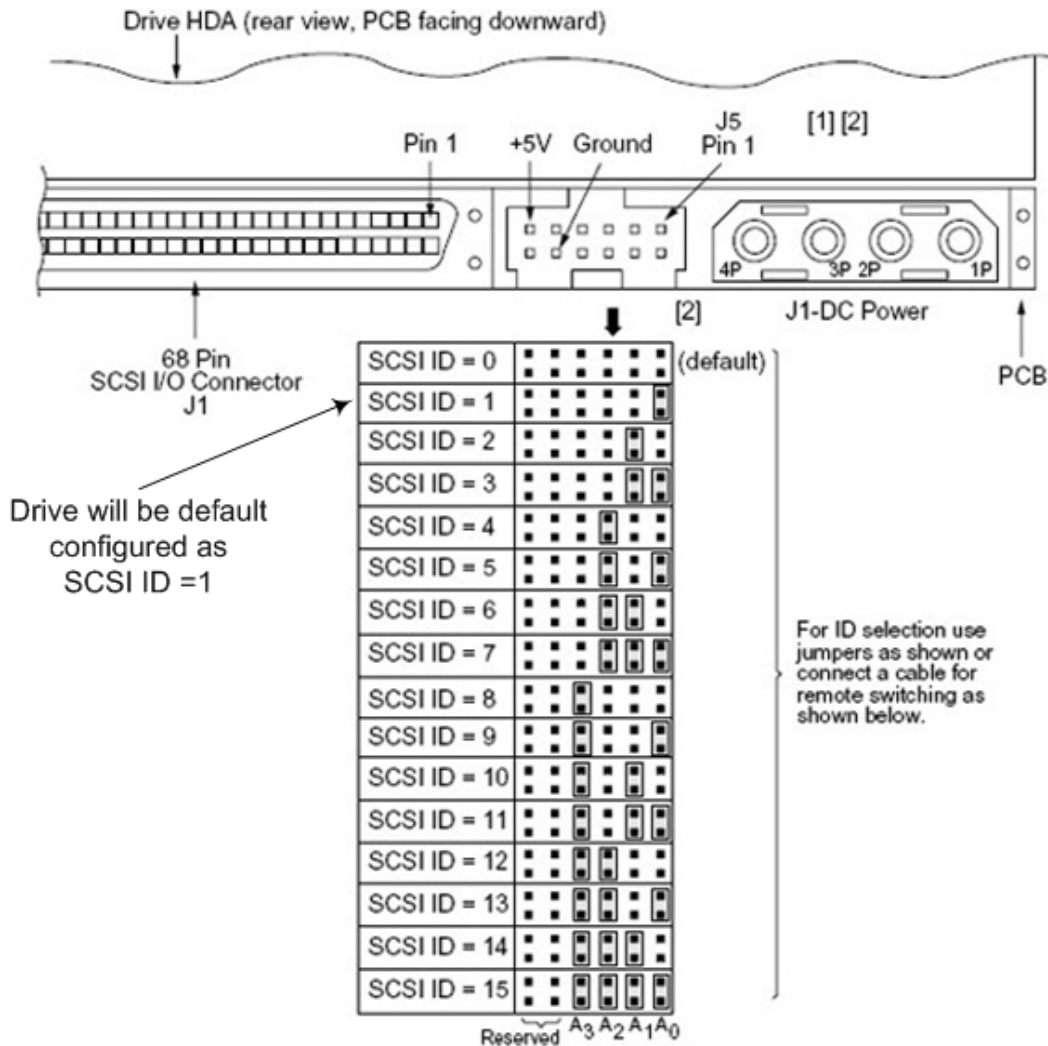

GE Healthcare MI&CT Class C eService Notes are copyright (c) protected to General Electric Company Inc with all rights reserved. Class C eService Notes issued by GE Healthcare are proprietary to GE Healthcare Systems. Disclosure to persons not employed by GE Healthcare or not in possession of a Advanced Service Agreement with GE Healthcare is prohibited.

Duplication and distribution of any class of eService Notes, in any form by entities other than GE Healthcare is strictly prohibited.

2225004-9 & -10 Hard Drives

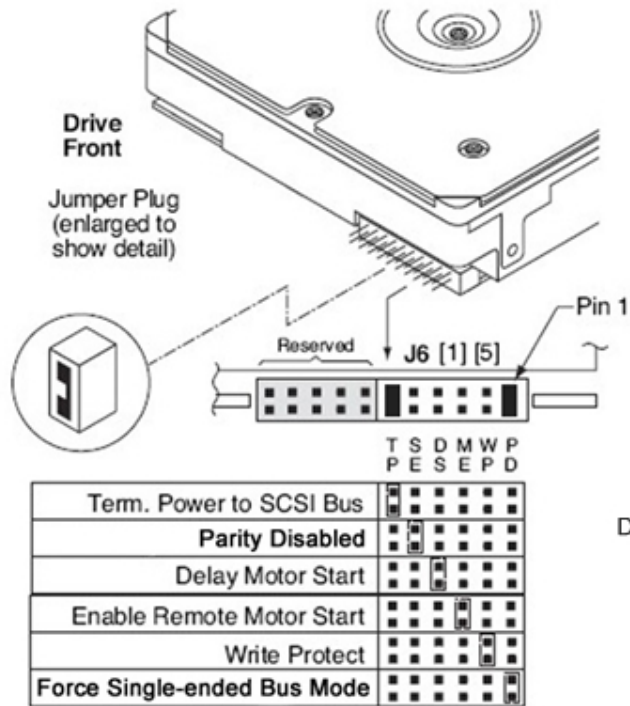

Drive will be default configured with
Force Single-ended Bus Mode
and
Term. Power to SCSI Bus
Jumpers Installed

GE Healthcare MI&CT Class C eService Notes are copyright (c) protected to General Electric Company Inc with all rights reserved. Class C eService Notes issued by GE Healthcare are proprietary to GE Healthcare Systems. Disclosure to persons not employed by GE Healthcare or not in possession of a Advanced Service Agreement with GE Healthcare is prohibited.
Duplication and distribution of any class of eService Notes, in any form by entities other than GE Healthcare is strictly prohibited.

Affected Service Publication:

LS-1X Jumpers 2282386-100, LS-1X IPL 2211225-100, LS-2X Jumpers 2243315-100, LS-2X IPL 2243316-100, and LS-3X H3, GOC1 and GOC2 Consoles IPLs and Jumper

PQR/SPR #: 13135055, 13103329, and 13103330

Note: FEMC "Parts Comments" have been updated to reflect latest part numbers.

Michael L. Bauer
Console Subsystem SMDE
MI&CT Service Engineering