



2-1 Cable map (Type A Detachable Table)

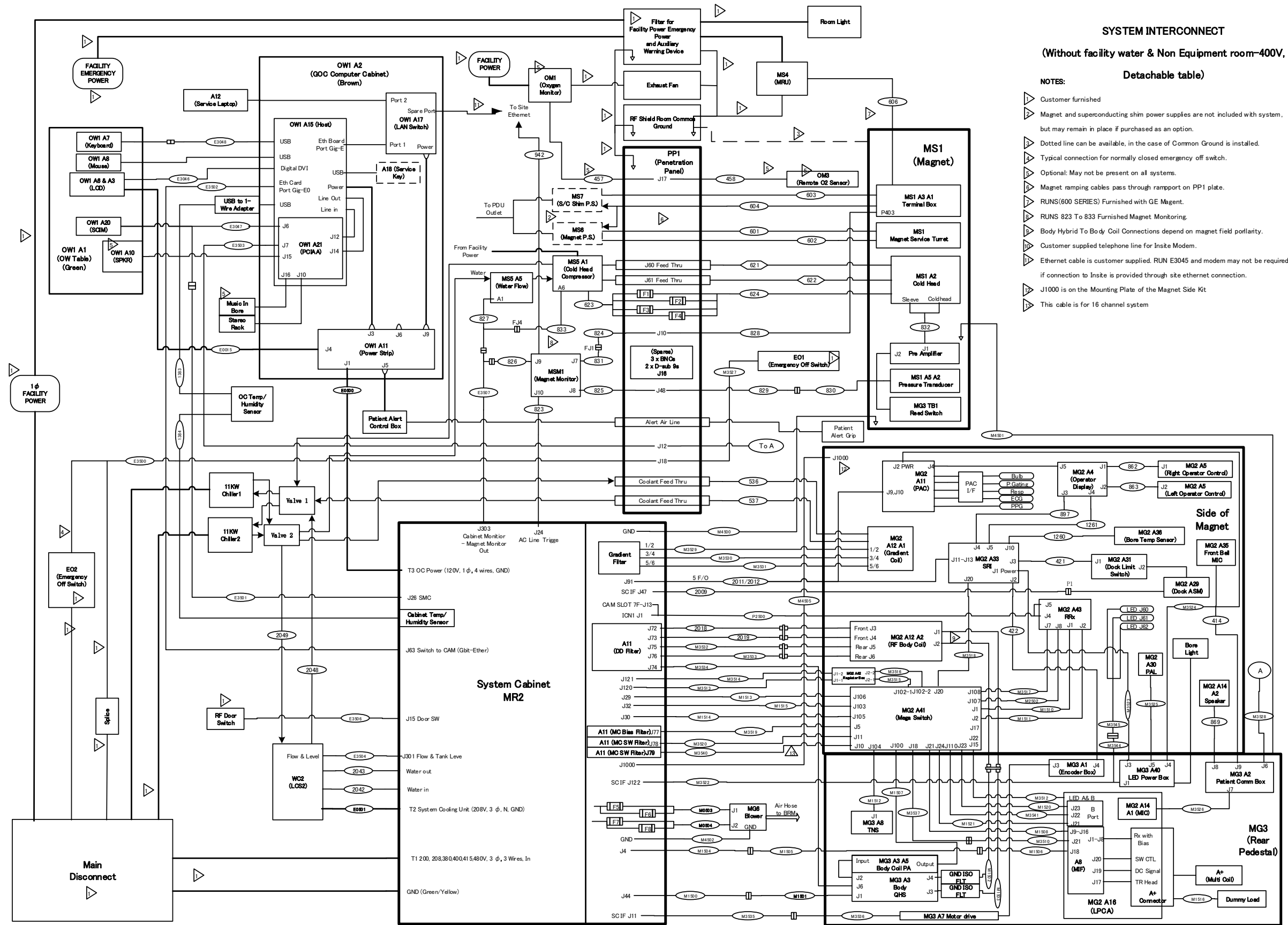

2-2 Cable map (Type B Detachable Table)

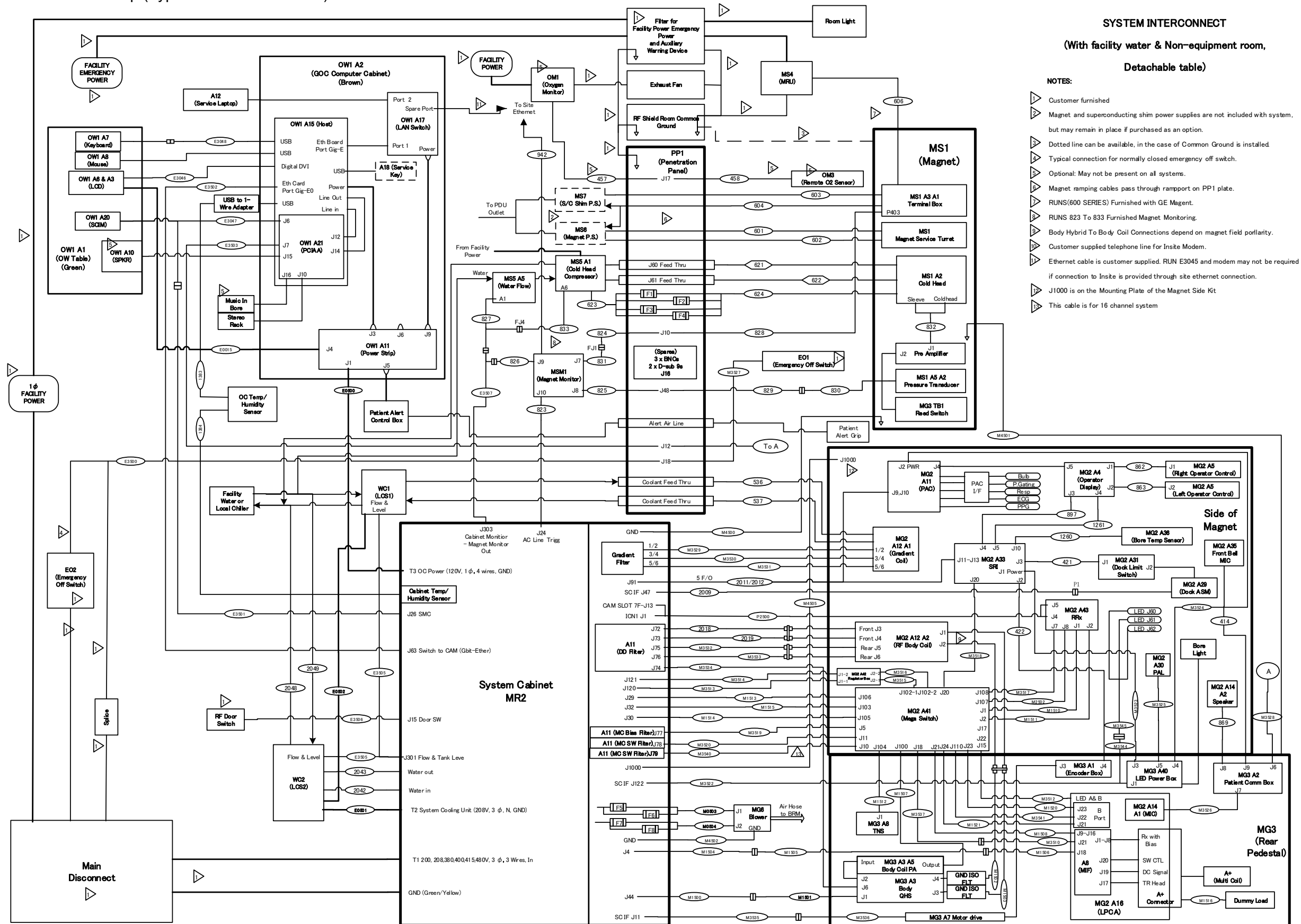


SYSTEM INTERCONNECT

(With facility water & Non-equipment room,
Detachable table)

NOTES:

- 1 Customer furnished
- 2 Magnet and superconducting shim power supplies are not included with system, but may remain in place if purchased as an option.
- 3 Dotted line can be available, in the case of Common Ground is installed.
- 4 Typical connection for normally closed emergency off switch.
- 5 Optional: May not be present on all systems.
- 6 Magnet ramping cables pass through ramping on PP1 plate.
- 7 RUNS(600 SERIES) Furnished with GE Magnet.
- 8 RUNS 823 To 833 Furnished Magnet Monitoring.
- 9 Body Hybrid To Body Coil Connections depend on magnet field polarity.
- 10 Customer supplied telephone line for InSite Modem.
- 11 Ethernet cable is customer supplied. RUN E3045 and modem may not be required if connection to InSite is provided through site ethernet connection.
- 12 J1000 is on the Mounting Plate of the Magnet Side Kit
- 13 This cable is for 16 channel system

2-3 Cable map (Type C Detachable Table)

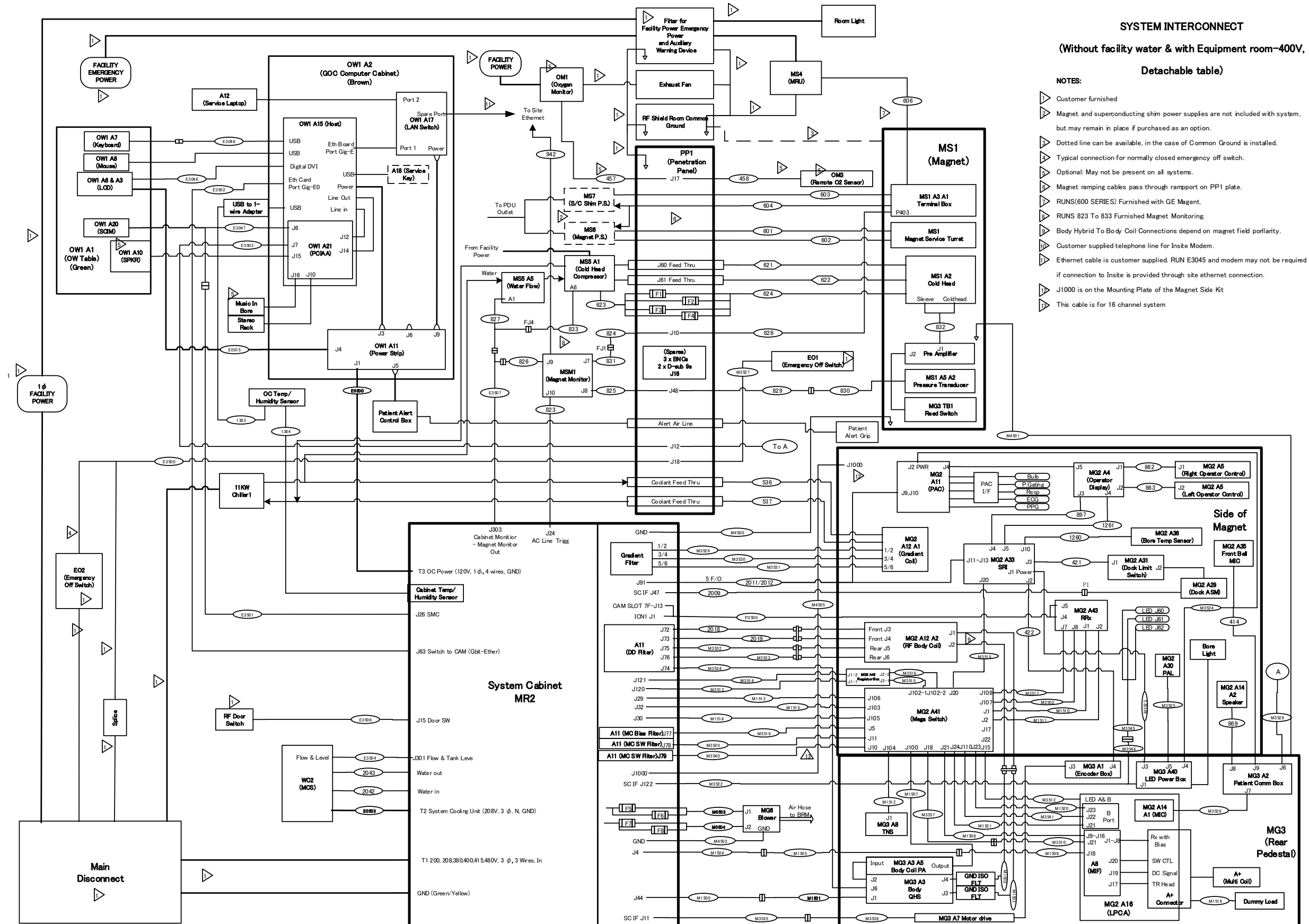


SYSTEM INTERCONNECT

(Without facility water & with Equipment room-400V,
Detachable table)


NOTES:

- ▶ Customer furnished
- ▶ Magnet and superconducting shim power supplies are not included with system, but may remain in place if purchased as an option.
- ▶ Dotted line can be available, in the case of Common Ground is installed.
- ▶ Typical connection for normally closed emergency off switch.
- ▶ Optional: May not be present on all systems.
- ▶ Magnet ramping cables pass through rampport on PP1 plate.
- ▶ RUNS(600 SERIES) Furnished with GE Magnet.
- ▶ RUNS 823 To 833 Furnished Magnet Monitoring.
- ▶ Body Hybrid To Body Coil Connections depend on magnet field polarity.
- ▶ Customer supplied telephone line for Insite Modem.
- ▶ Ethernet cable is customer supplied. RUN E3045 and modem may not be required if connection to Insite is provided through site ethernet connection.
- ▶ J1000 is on the Mounting Plate of the Magnet Side Kit
- ▶ This cable is for 16 channel system


2-4 Cable map (Type D Detachable Table)

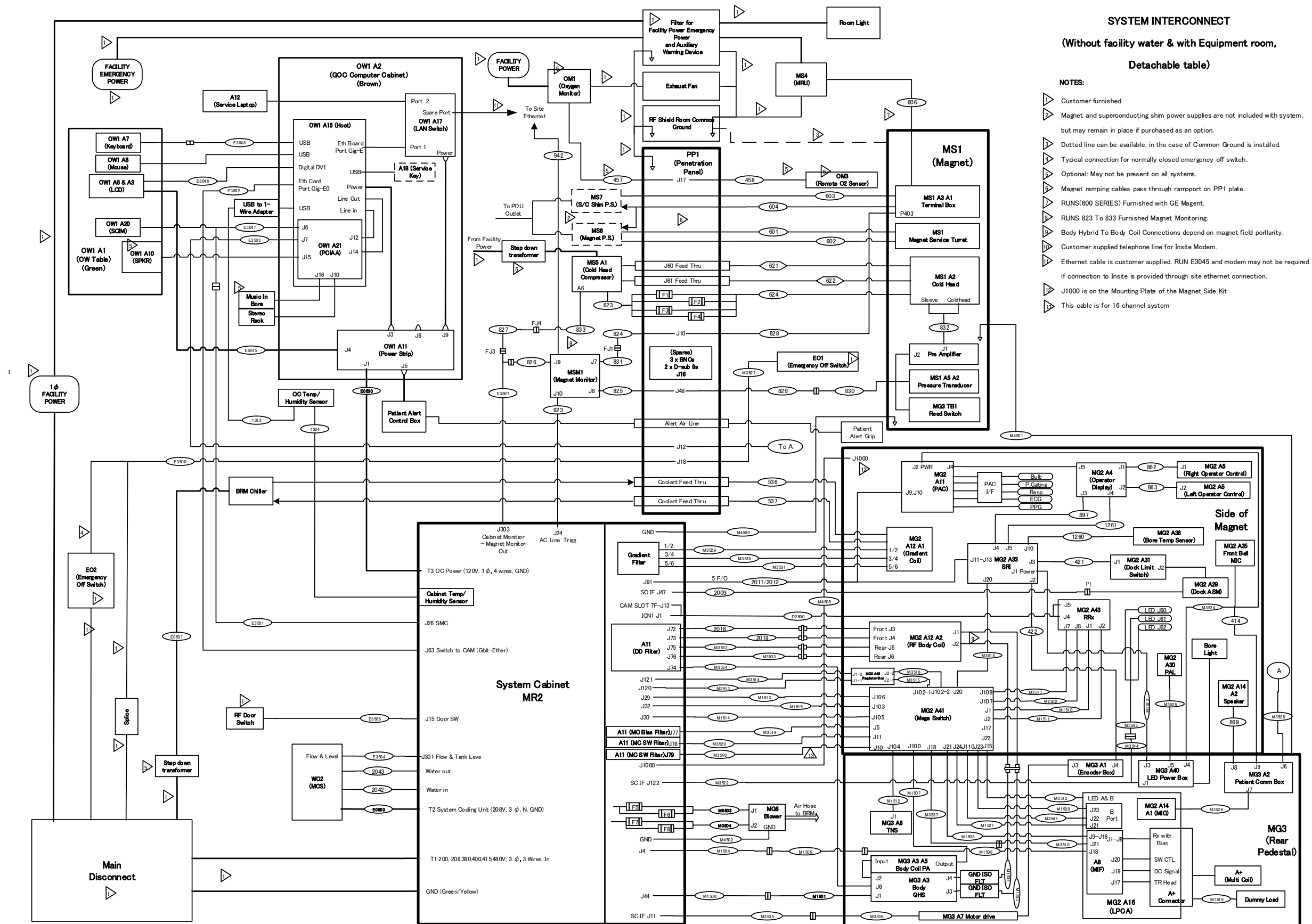


2-5 Cable map (Type E Detachable Table)

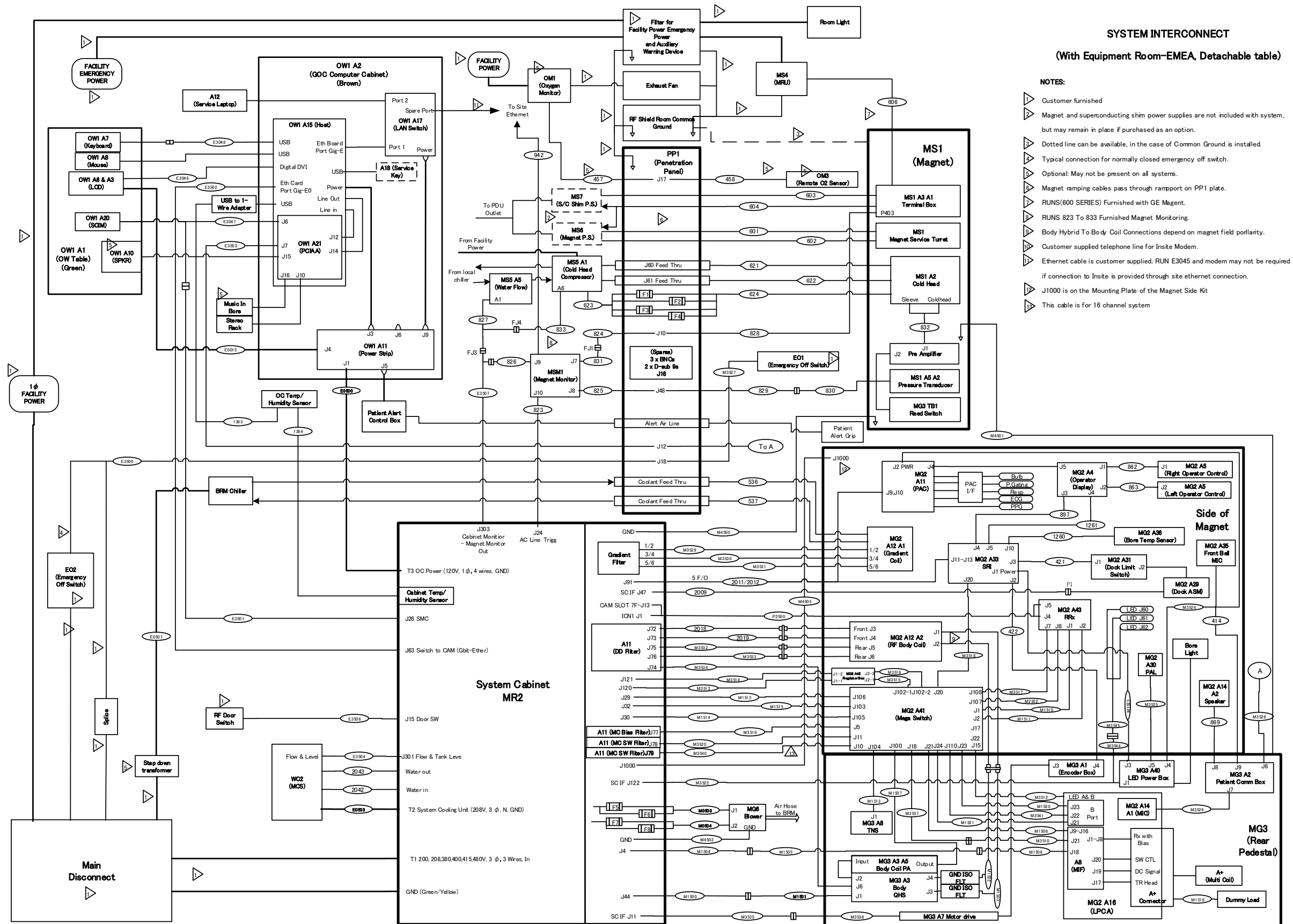

2-6 Cable map (Type A Fixed Table)

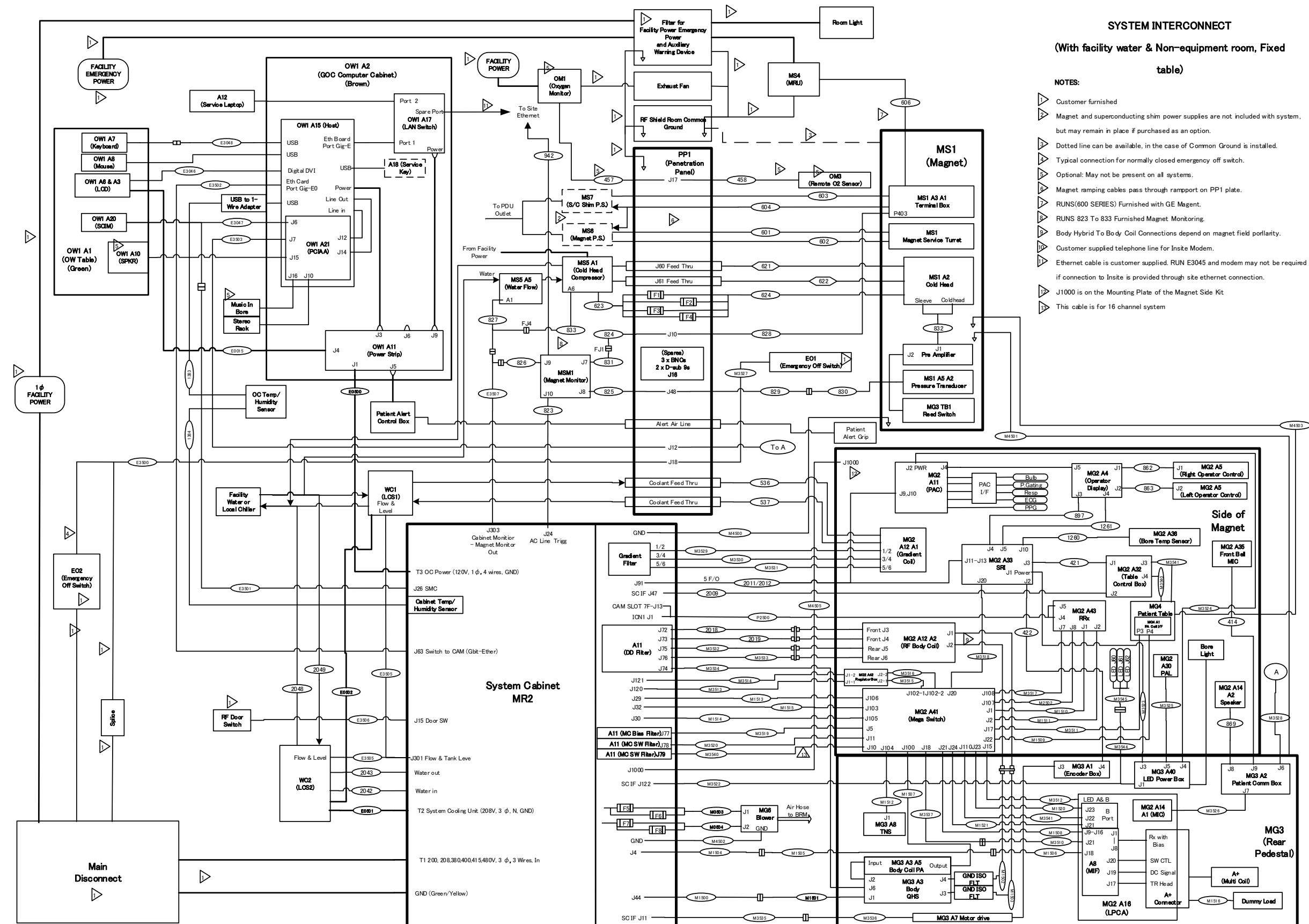


2-7 Cable map (Type B d Table)

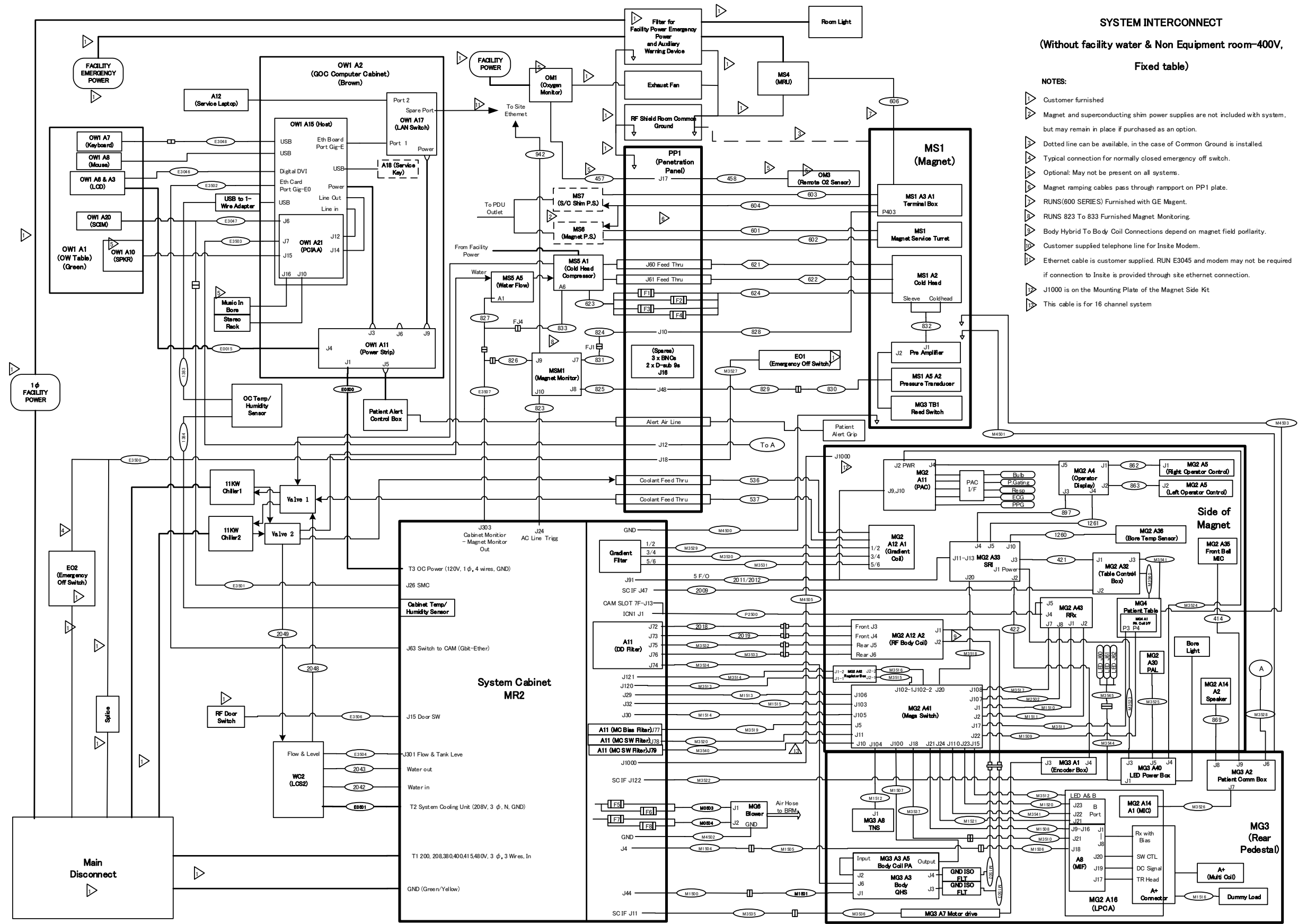


SYSTEM INTERCONNECT
(Without facility water & Non Equipment room-400V,
Fixed table)

- NOTES:
- Customer furnished
 - Magnet and superconducting shim power supplies are not included with system, but may remain in place if purchased as an option.
 - Dotted line can be available, in the case of Common Ground is installed
 - Typical connection for normally closed emergency off switch.
 - Optional: May not be present on all systems.
 - Magnet ramping cables pass through rampport on PP1 plate.
 - RUNS(600 SERIES) Furnished with GE Magnet.
 - RUNS 823 To 833 Furnished Magnet Monitoring.
 - Body Hybrid To Body Coil Connections depend on magnet field portlarity.
 - Customer supplied telephone line for Insite Modem.
 - Ethernet cable is customer supplied. RUN E3045 and modem may not be required if connection to Insite is provided through site ethernet connection.
 - J1000 is on the Mounting Plate of the Magnet Side Kit
 - This cable is for 16 channel system

2 8 Cable map (Type C Fixed Table)

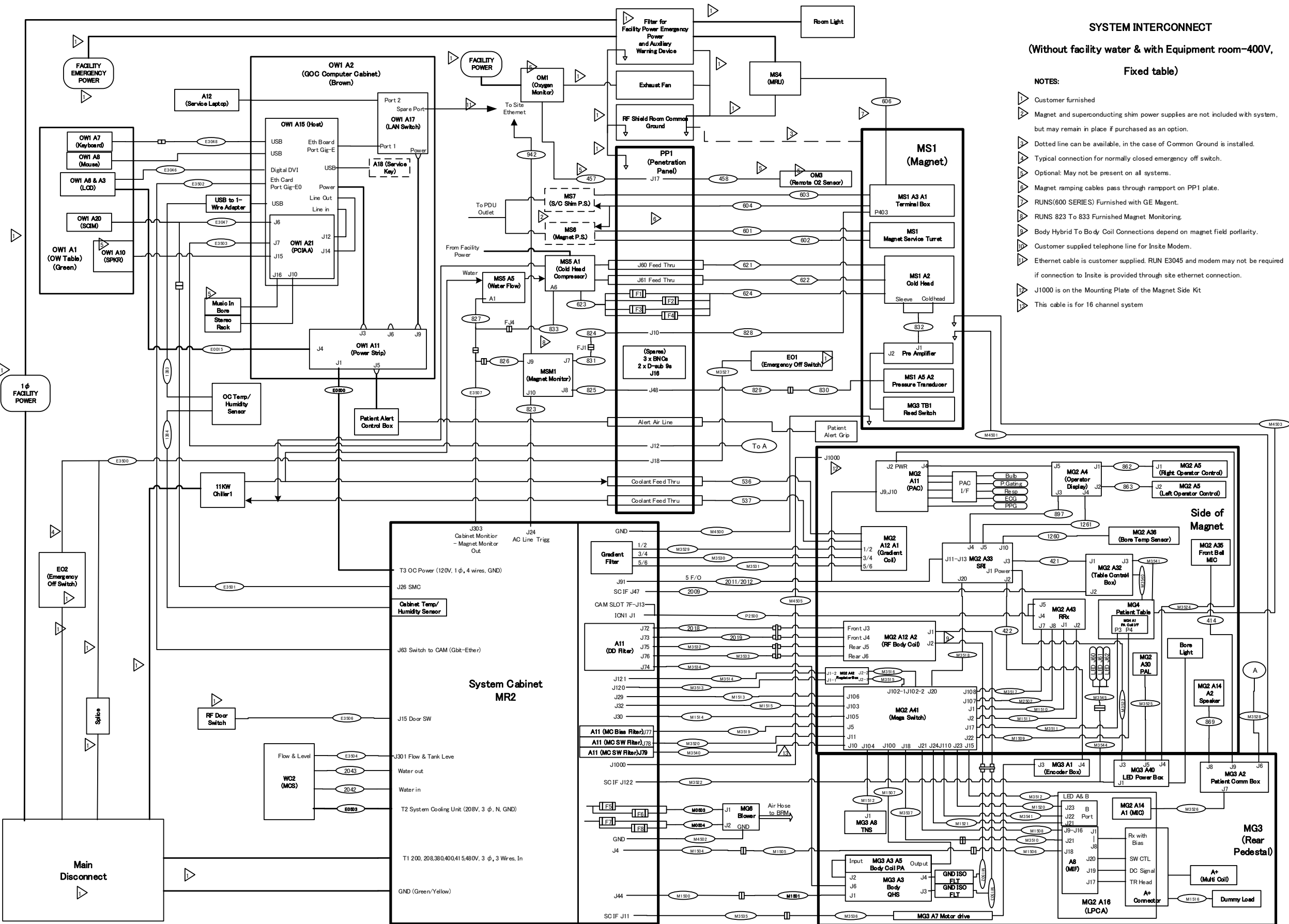

2-9 Cable map (Type D Fixed Table)

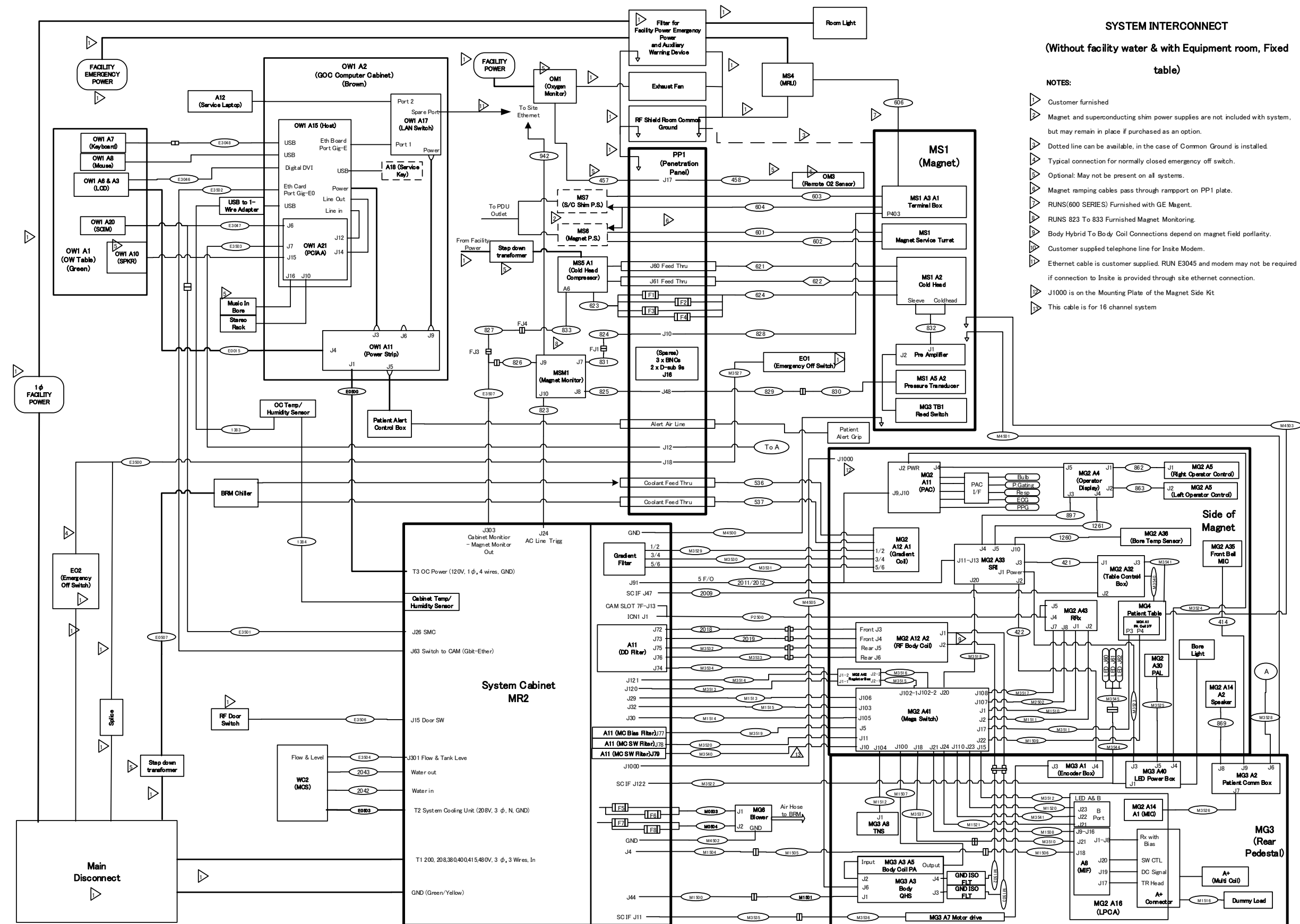

2- Cable map (Type Fixed Table)

