Technical Publications

Direction 5121980
Revision 1

Signa® ReFlex AP Option Installation for Excite HD

Copyright© 2004 by General Electric Company, Inc.
All Rights Reserved

Operating Documentation

DAMAGE IN TRANSPORTATION

All packages should be closely examined at time of delivery. If damage is apparent, have notation "**damage in shipment**" written on **all** copies of the freight or express bill **before** delivery is accepted or "signed for" by a General Electric representative or a hospital receiving agent. Whether noted or concealed, damage **MUST** be reported to the carrier **immediately** upon discovery, or in any event, within **14** days after receipt, and the contents and containers held for inspection by the carrier. A transportation company will not pay a claim for damage if an inspection is not requested within this **14** day period.

Immediately complete a "Damage Loss Claim Form", available via MS Exchange Mail, after the damage is found.

MS Exchange Path:

Outlook/Public Folder/All Public Folders/Medical Systems/!Global Initiatives/Information Management/Forms/Common Forms/DAMAGE LOSS CLAIM FORM.

Send the completed form to the email address listed in the form.

For more information about the Transportation Claim Procedure, access the GE Medical Systems Intranet and enter the following URL address (case sensitive):

<ftp://3.87.40.2/globepro/qualsys/Docs/190016MF.PDF>

Rev. 11/15/2000

Direction 2128126 - Language Policy For Service Documentation, Rev 0

WARNING

- THIS SERVICE MANUAL IS AVAILABLE IN ENGLISH ONLY.
- IF A CUSTOMER'S SERVICE PROVIDER REQUIRES A LANGUAGE OTHER THAN ENGLISH, IT IS THE CUSTOMER'S RESPONSIBILITY TO PROVIDE TRANSLATION SERVICES.
- DO NOT ATTEMPT TO SERVICE THE EQUIPMENT UNLESS THIS SERVICE MANUAL HAS BEEN CONSULTED AND IS UNDERSTOOD.
- FAILURE TO HEED THIS WARNING MAY RESULT IN INJURY TO THE SERVICE PROVIDER, OPERATOR OR PATIENT FROM ELECTRIC SHOCK, MECHANICAL OR OTHER HAZARDS.

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* This revision/letter corresponds to the indicated document's revision control system.

TABLE OF CONTENTS

1- INTRODUCTION	7
1-1 Signa 12.x Options Installed by this Manual	7
1-2 Contents	7
2- REFLEX AP INSTALLATION (NEW INSTALLATIONS).....	8
2-1 Power Off Procedure	8
2-2 Installation Reflex 200B	9
2-2-1 Set-up of Master Slave boards	9
2-2-2 Installation of Option	10
2-3 Installation of Reflex 400B or 800	11
2-3-1 Reflex 400B installation	11
2-3-2 Reflex 800 installation.....	11
2-4 Power On and TPS Reset	12
2-5 Functional Check	13

This document is only compatible for systems with Excite HD (or 12x) software.

1- INTRODUCTION

This document covers new installations and upgrades of the Mercury AP boards with Excite HD.

Total installation time for the hardware is up to one hour for one installer.

1-1 Signa 12.x Options Installed by this Manual

- M3000RL Excite ReFlex 200B
- M3000RM Excite ReFlex 400B
- M3000RN Excite ReFlex 800

1-2 Contents

M3000RL (ReFlex200B):

- **Two** motherboards 2325947-11
- Rating Plate
- PIM Module

M3000RM (ReFlex400B):

- Minerva II Master board 2325947-12
- Minerva II Slave board 2325947-13
- Rating Plate
- PIM Module

M3000RN (ReFlex800):

- Minerva II Master board 2325947-15
- Minerva II Slave board 2325947-16
- Rating Plate
- PIM Module


2- REFLEX AP INSTALLATION (NEW INSTALLATIONS)

EQUIPMENT DAMAGE POSSIBILITY. Electro Static Discharge (ESD) precautions must be used when removing the AP assembly or replacing any modules on the AP.

2-1 Power Off Procedure

1. Remove the System Cabinet front cover and open rear door.
2. Turn off power switch at rear of MGD chassis. See Illustration 2-1.

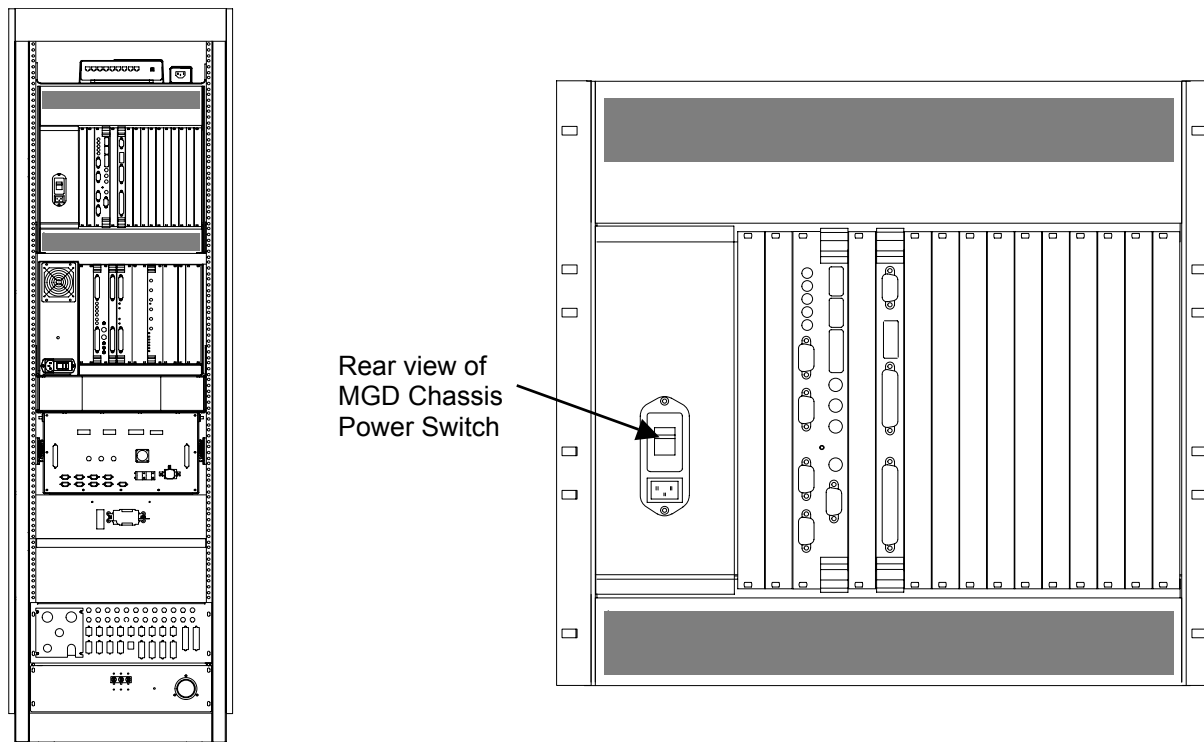
SYSTEM CABINET POWER OFF LOCATIONS
ILLUSTRATION 2-1

2-2 Installation of Reflex 200B

2-2-1 Set-up of Master Slave boards

Rotary Switch SW2 is used to determine Master or Slave boards. See Illustration 2-2 for the location of SW2.


The **Master** board must have SW2 set to **0** should be installed in slot 2.
The **Slave** board must have SW2 set to **2** should be installed in slot 3.

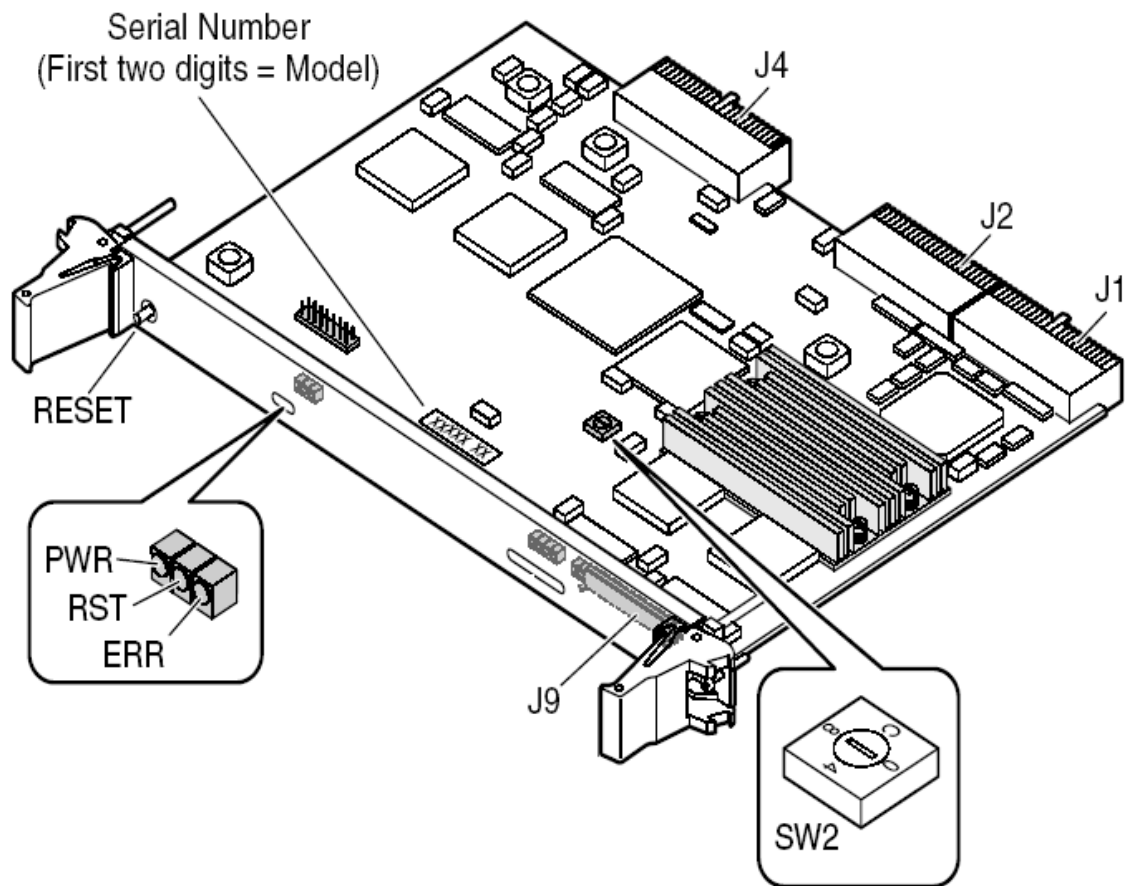
MASTER SLAVE SELECTION
ILLUSTRATION 2-2

2-2-2 Installation of Option

1. Install Master board in slot 2, then install Slave board in slot 3.
2. Proceed to section 2-4, installation of PIM module.

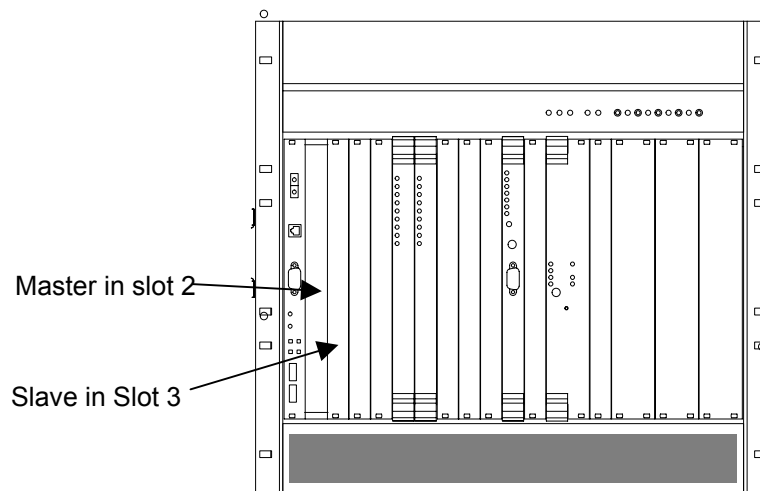BOARD LOCATIONS
ILLUSTRATION 2-3

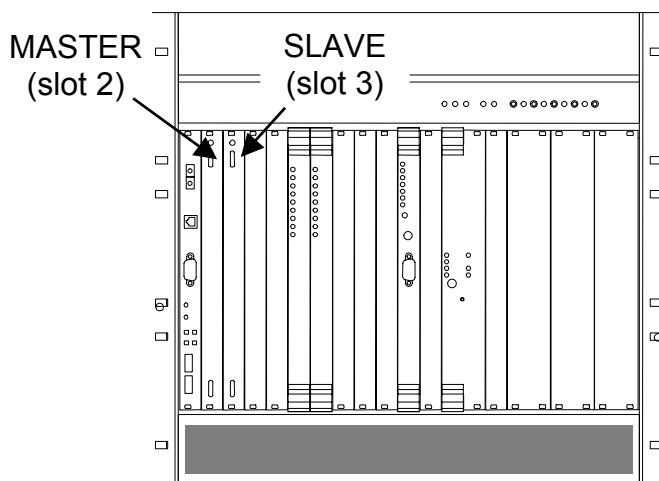NOTE:
The SLAVE board must always be to the right of the MASTER board.

INSTALLATION OF REFLEX200B IN MGD CHASSIS
ILLUSTRATION 2-4

2-3 Installation of ReFlex400B or ReFlex800

2-3-1 ReFlex 400B installation

The ReFlex400B is made up of two boards with two daughter-cards on each board. There are no adjustments or jumper settings are necessary for the Reflex400B. Also there is no front raceway needed for Reflex400B.

1. Insert part number **2325947-12** (label 900-07072) into slot 2. See Illustration 2-3.
2. Insert part number **2325947-13** (label 900-07071) into slot 3.
3. Fasten in place with screws at top and bottom of front plate.
4. Go to PIM Module Installation, section 2-4.

2-3-2 ReFlex 800 installation

The ReFlex800 is made up of two boards with two daughter-cards on each board. There are no adjustments or jumper settings are necessary for the Reflex800. Also there is no front raceway needed for Reflex800.

1. Insert part number **2325947-15** (label 900-07075) into slot 2. See Illustration 2-3.
2. Insert part number **2325947-16** (label 900-07078) into slot 3.
3. Fasten in place with screws at top and bottom of front plate.
4. Go to PIM Module Installation, section 2-4.

Note

The backplane of the Excite II cabinet contains the necessary jumpers between slots 2 and 3. Do not place the raceway adaptor or anything between the two boards on the front of the MGD chassis.

2-4 Installation of PIM module

Note

Some cabinets may already have PIM module installed. If the PIM is installed, proceed to step 2.5.

Note

Use caution when inserting the PIM module. The backplane pins can be damaged if excessive force is used or if the PIM is not straight when being plugged in.

1. Remove four faceplates from the rear of the MGD chassis. See Illustration 2-5.
2. Install PIM module onto the MGD backplane between slots 2 and 3. See Illustration 2-5 and 2-6. The PIM module allows communication between the Master and Slave boards.

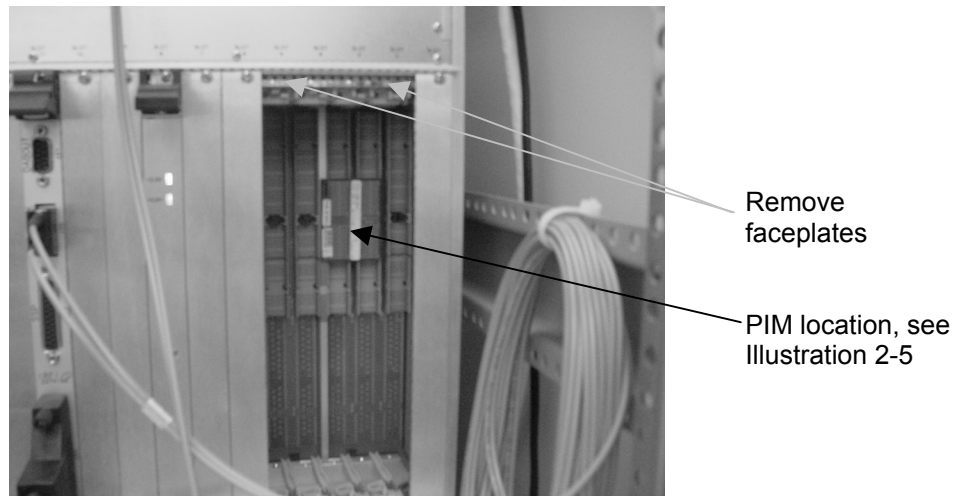
Illustration 2-5
PIM Installation

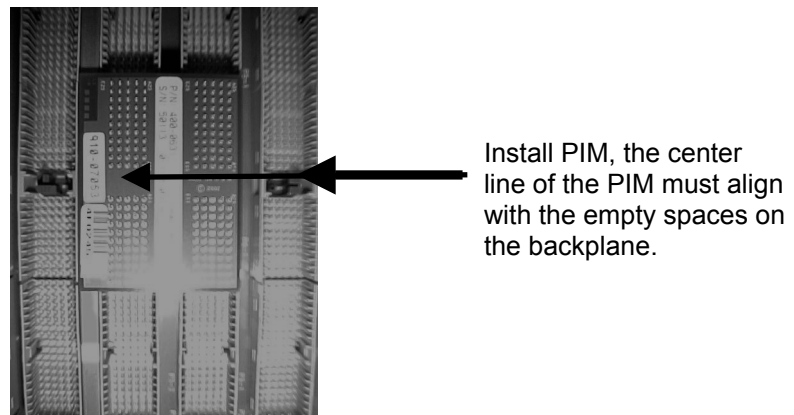
Illustration 2-6
PIM Installation

3. Re-install the four faceplates back onto the MGD chassis.

2-5 Power On and TPS Reset

After installing the AP board(s), power on the MGD chassis, then do a **[TPS Reset]**.

There is no need to modify the configuration file because the system software will auto-configure the new board(s).

2-6 Functional Check

Run AP diags by doing the following.

1. Open the Common Service desktop or CSD.
2. Select **Manual Diags**, then select **AP functional check**.
3. Select **Run**. You should then see which AP boards were recognized by the system.
See Illustration 2-7.

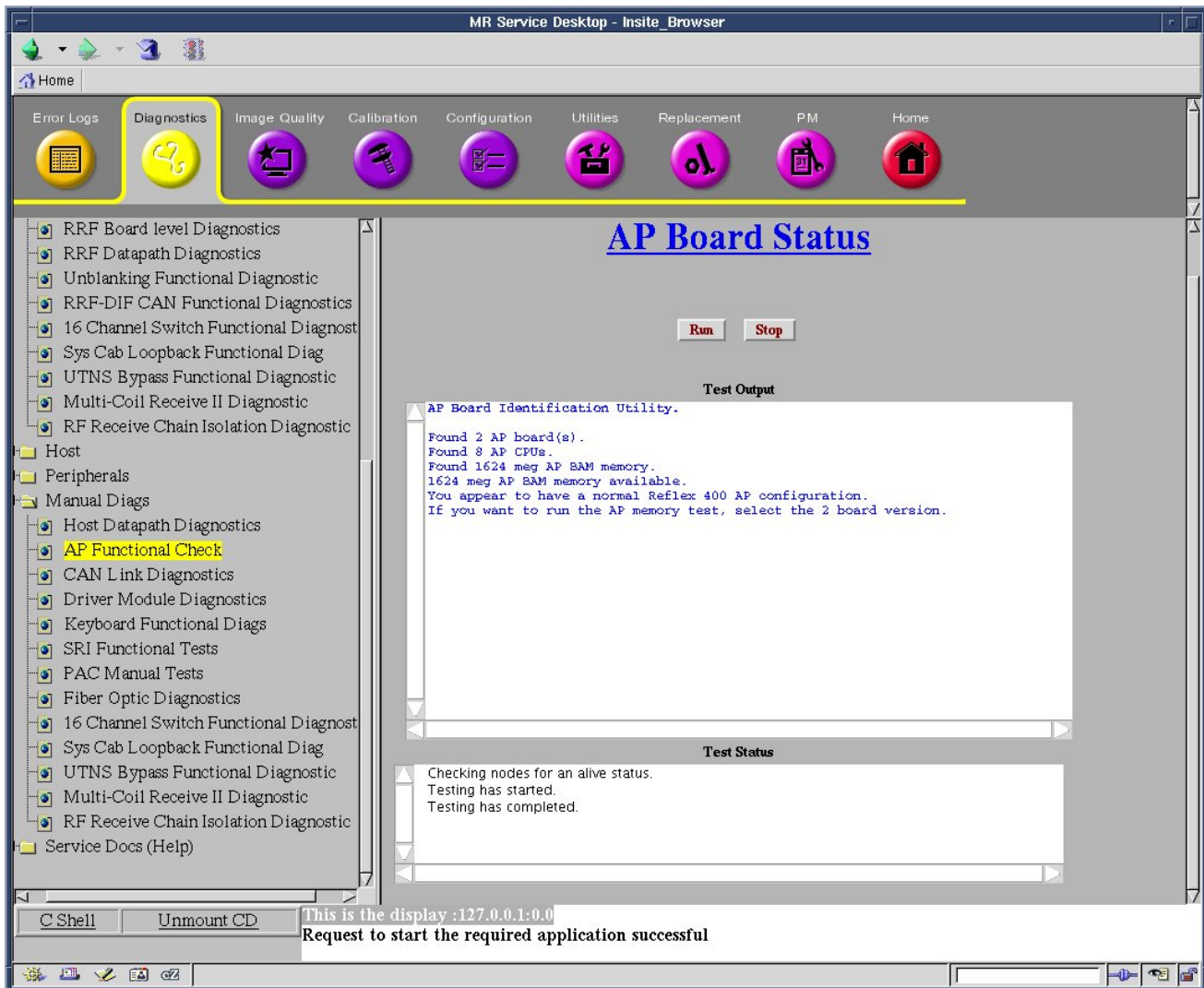
Illustration 2-7