GE Medical Systems

Technical Publications

**Direction 2209382
Revision 6**

Signa[®] Software Option Installation

**Copyright © 1998, 1999, 2000, 2001, 2003 by General Electric Company, Inc.
All Rights Reserved**

Operating Documentation

DAMAGE IN TRANSPORTATION

All packages should be closely examined at time of delivery. If damage is apparent, have notation "**damage in shipment**" written on **all** copies of the freight or express bill **before** delivery is accepted or "signed for" by a General Electric representative or a hospital receiving agent. Whether noted or concealed, damage **MUST** be reported to the carrier **immediately** upon discovery, or in any event, within **14** days after receipt, and the contents and containers held for inspection by the carrier. A transportation company will not pay a claim for damage if an inspection is not requested within this **14** day period.

Immediately complete a "Damage Loss Claim Form", available via MS Exchange Mail, after the damage is found.

MS Exchange Path:

Outlook/Public Folder/All Public Folders/Medical Systems/!Global Initiatives/Information Management/Forms/Common Forms/DAMAGE LOSS CLAIM FORM.

Send the completed form to the email address listed in the form.

For more information about the Transportation Claim Procedure, access the GE Medical Systems Intranet and enter the following URL address (case sensitive):

<ftp://3.87.40.2/globepro/qualsys/Docs/190016MF.PDF>

Rev. 11/15/2000

Direction 2128126 - Language Policy For Service Documentation, Rev 0

W A R N I N G

- THIS SERVICE MANUAL IS AVAILABLE IN ENGLISH ONLY.
- IF A CUSTOMER'S SERVICE PROVIDER REQUIRES A LANGUAGE OTHER THAN ENGLISH, IT IS THE CUSTOMER'S RESPONSIBILITY TO PROVIDE TRANSLATION SERVICES.
- DO NOT ATTEMPT TO SERVICE THE EQUIPMENT UNLESS THIS SERVICE MANUAL HAS BEEN CONSULTED AND IS UNDERSTOOD.
- FAILURE TO HEED THIS WARNING MAY RESULT IN INJURY TO THE SERVICE PROVIDER, OPERATOR OR PATIENT FROM ELECTRIC SHOCK, MECHANICAL OR OTHER HAZARDS.

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REVISION HISTORY

<u>REV</u>	<u>DATE</u>	<u>PRIMARY REASON FOR CHANGE</u>
0.....	Feb 11, 1998	Initial release.
1.....	June 4, 1998	Updated for release 8.2
2.....	June 15, 1999	Updated for release 8.2.5
3.....	May, 2000	Conversion Interleaf to Word. Added Catalog Availability Info for 8.2
4.....	September, 2000	Updated Catalog Information to include Release 8.3 Information.
5.....	August, 2001	Updated Catalog information to include ASP, 9.0 and OpenSpeed. Revised Document name to be more generic. Added a Glossary.
6.....	August 4, 2003	Added information for 11x Floppy load. Removed part numbers, they do not belong in the installation document

<u>PAGE</u>	<u>REV</u>	<u>PAGE</u>	<u>REV</u>	<u>PAGE</u>	<u>REV</u>	<u>PAGE</u>	<u>REV</u>
Title Page.....	6						
2 of 8 to 8 of 8	6						

* This revision/letter corresponds to the indicated document's revision control system.

TABLE OF CONTENTS

1 OVERVIEW..... 7

2 OPTION SOFTWARE INSTALLATION 7

3 REBOOT COMPUTER 8

1 OVERVIEW

This procedure describes how to install software options on the Signa scanners 8.3, ASP, 9.x, 10x, 11x or later systems using the MOD or Floppy drive. With the introduction of 11x, Software Option Keys are provided via floppy. The Software Option package contains the appropriate software option MOD or Floppy.

2 OPTION SOFTWARE INSTALLATION

1. Media Installation
 - a. Insert the MOD disk in the SGI computer MOD drive for 10x of previous revisions.
 - b. Insert the Floppy into the Linux PC for 11x sites.
2. Select the **[Common Service Desktop]**.
3. Select **[Configuration]** from the Service Desktop Manager.
4. Select **[Guided Install]**
5. Select FE Mode if loading to Signa software revision ASP or higher.
6. A command window will open with "**Type Root Password**" displayed. Type the systems root password and then press the **<Enter>** key. "**Operator**" is the Signa default root password. If for your site, the root password has been changed by the customer you will have to obtain the correct root password from the customer.
7. After the Guided Install window opens (10 sec for Linux PC, 30 sec. for Octane; 60 sec. for Indigo), select **Option Keys** tab. See Illustration 2-1.

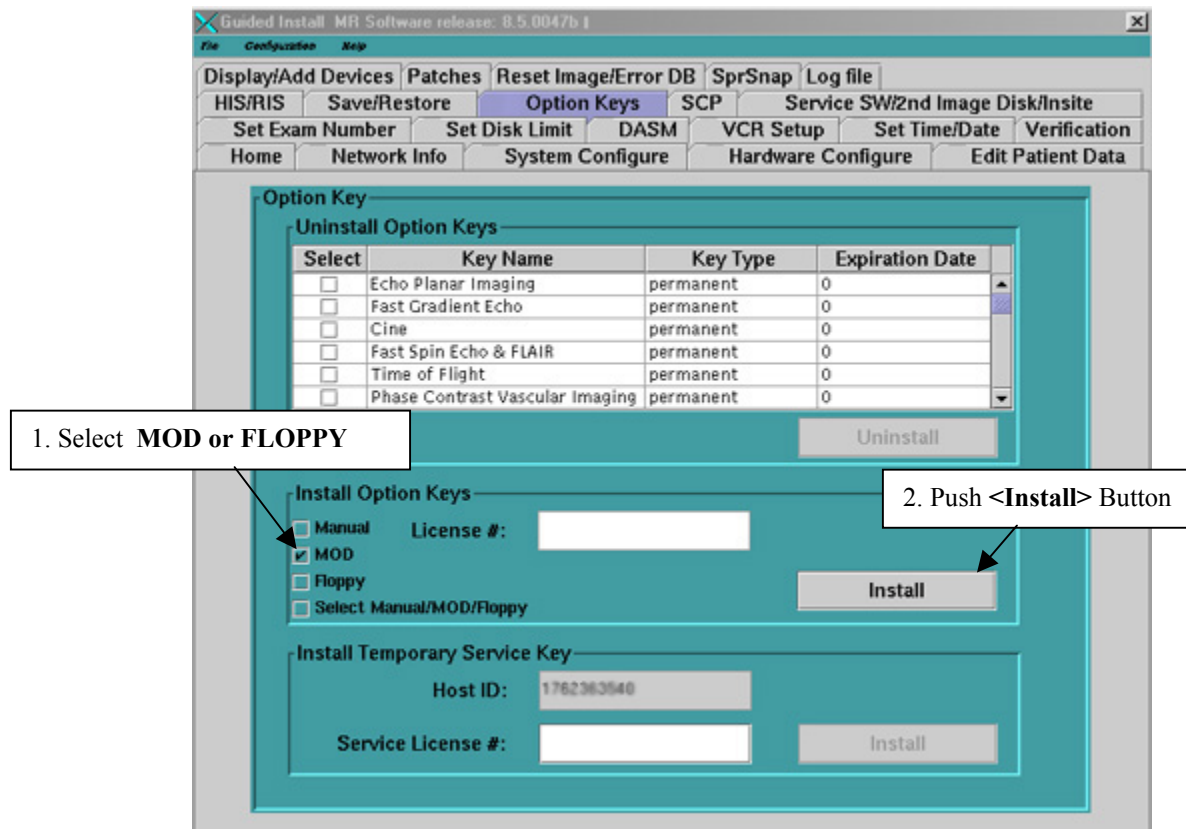
ILLUSTRATION 2-1
GUIDED INSTALL WINDOW

8. Select **MOD** or the **Floppy** box at Install Option Key portion of window, then select **Install**.
9. The Option Key will begin loading. Load times for option keys are very short. Typically less than 1 minute. When the Option Key is finished loading a popup will appear requesting you to insert the next Option Key (MOD or floppy) or press **<Quit>**. See Illustration 2-2. If other options are to be loaded at this time, remove the currently loaded Option Key by pressing the eject button on the Option Key drive and insert the next Option Key to be loaded. Wait 15 seconds for loading to begin. Continue this process until all options are loaded.

ILLUSTRATION 2-2
INSTRUCTIONAL POP UP FOR OPTION KEY LOADS

10. After all options are loaded, press the **<Quit>** button to quit the option load process. Select **File ->Quit** from the main menu bar to exit the Guided Install tool.

3 REBOOT COMPUTER

1. Reboot the computer.
2. The installation is now complete.