



Downscan Converter Kit Configuration Direction

This direction includes connection diagrams for video
downscanner configurations in the following kits:

S9200LT

Advantage Windows Workstation 2.0/3.1	pg. 4
MR Contour/Profile 3 (color video)	pg. 5
MR Cardio/Neuro Vascular	pg. 5
MR Lightning LX	pg. 6
Xray GEMNet	pg. 6, 7
MR PC Linux Host	pg. 8

S9200LW

MR Contour/Profile (Monochrome Video)	pg. 8
---------------------------------------	-------

S9200LS

Xray Prestige/Advantix R&F, Vascular, Cardiac (DRS/DLX)	pg. 9
---	-------


S9200LY

MR Lightning LX PV Replacement	pg. 9
--------------------------------	-------

Revision History

<u>Date/Rev</u>	<u>Changes</u>	<u>Author</u>	<u>Approved</u>
03/04/1999 Rev. 1.9	Pg 4: Swapped DIP settings for AW3.1 MS8154 & MS8153 cables.	J. Rutschky	J. Rutschky
03/22/1999 Rev. 2.0	Pg 8: Corrected typo error, : MS8165CA to MS8156CA. Added HP B180L cable configuration for E7200BJ.	J. Rutschky	J. Rutschky
05/12/0999 Rev. 3.0	All pages changed to reflect the consolidation of ten separate catalog items into three. The Sony printer (UP7300) will be configured as a separate optional kit for use with all GEMNet configurations.	J. Rutschky	J. Rutschky
09/23/1999 Rev. 3.1	E8003FG and E8003FH typographical errors were changed to E8803FG and E8803FH respectively.	J. Rutschky	J. Rutschky
08/08/2001 Rev. 4.0	E8803SA, E8803SB for MR Lightning LX added to replace the SGI Octane "Personal Video" option that is no longer available on SGI Octane II. Configuration diagram added on page 9.	J. Rutschky	J. Rutschky
05/13/2003 Rev. 4.1	Changed kit E-cat numbers to S-cat numbers to reflect new structuring. E8803FG changed to S9200LT E8803FH changed to S9200LW E8803EA changed to S9200LS E8803SA/SB changed to S9200LY	J. Rutschky	J. Rutschky
06/02/2003	New HP PC Linux Host drawing added to page 8 for the S9200LT configuration.	J. Rutschky	J. Rutschky
06132003	Changed two references for CB4012MR on page 3 to CB2012MR because the former was an incorrect Altinex part number.	J. Rutschky	J. Rutschky


Common to all configurations: Downscanner, VCR and Monitor

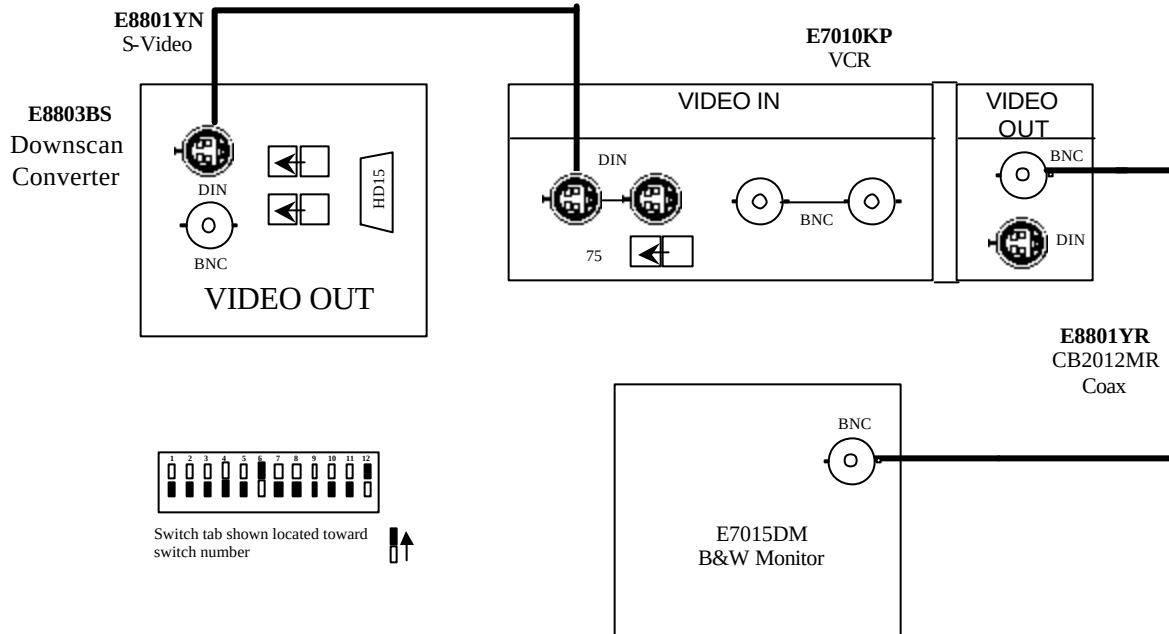


Parts List

<u>CAT #</u>	<u>Vendor #</u>		<u>Description</u>
E8803BS	DSC-1024		Downscan Converter
E7010KP	SVO-9500		Video Cassette Recorder
E7015DM	PVM-97		9" B&W Video Monitor
E7014PC	UPZ-7000I	Cable	Sony Printer-HD15M
E8801YA	DA1907SX		Video Dist. Amp.HD15
E8801YB	DA1907GBW		Video Dist. Amp.HD15
E8801YP	DA1804NT		Video Dist. Amp.BNC
E8801YC	MS8126CA	Cable	HD15M-HD15M 3'
E8801YD	MS8156CA	Cable	HD15M-HD15M 6'
E8801YE	MS8166CA	Cable	HD15M-HD15M 12'
E8801YF	MS8123CA	Cable	13W3F-HD15M 3'
E8801YG	MS8153CA	Cable	13W3F-HD15M 6'
E8801YH	MS8163CA	Cable	13W3F-HD15M 12'
E8801YJ	MS8124CA	Cable	13W3M-HD15M 3'
E8801YK	MS8154CA	Cable	13W3M-HD15M 6'
E8801YL	MS8164CA	Cable	13W3M-HD15M 12'
E8801YM	MS8102CA	Cable	HD15M-5BNM 6'
E8801YN	MS8160MG	Cable	S-Video (Y/C) 12'
E8801YR	CB2012MR	Cable	75 ohm BNC-BNC 12'
E7014PC	UPZ-7000I	Cable	High Scan Printer Video 3'
E8801Y	MS8146MG	Terminator BNC Barrel	HD15M-HD15F Supply Locally


S9200LT Advantage Workstation 2.0

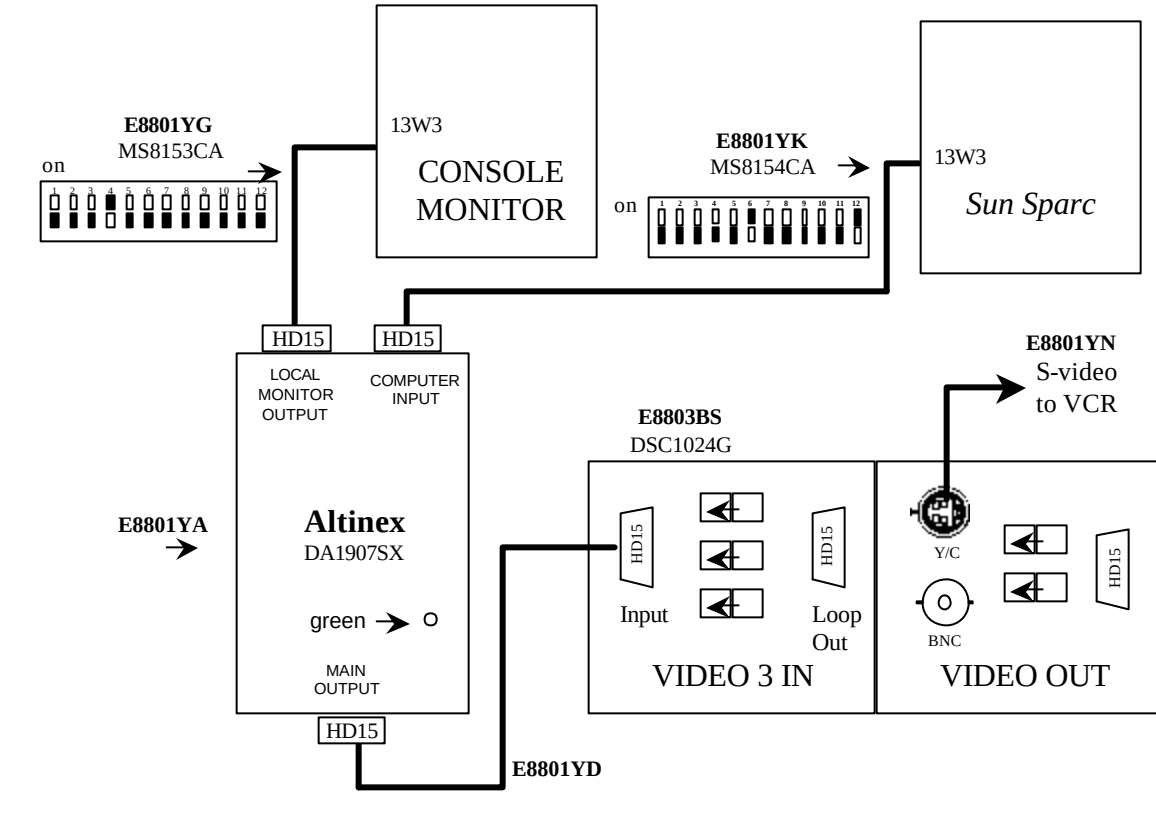


S9200LT Advantage Workstation 3.1

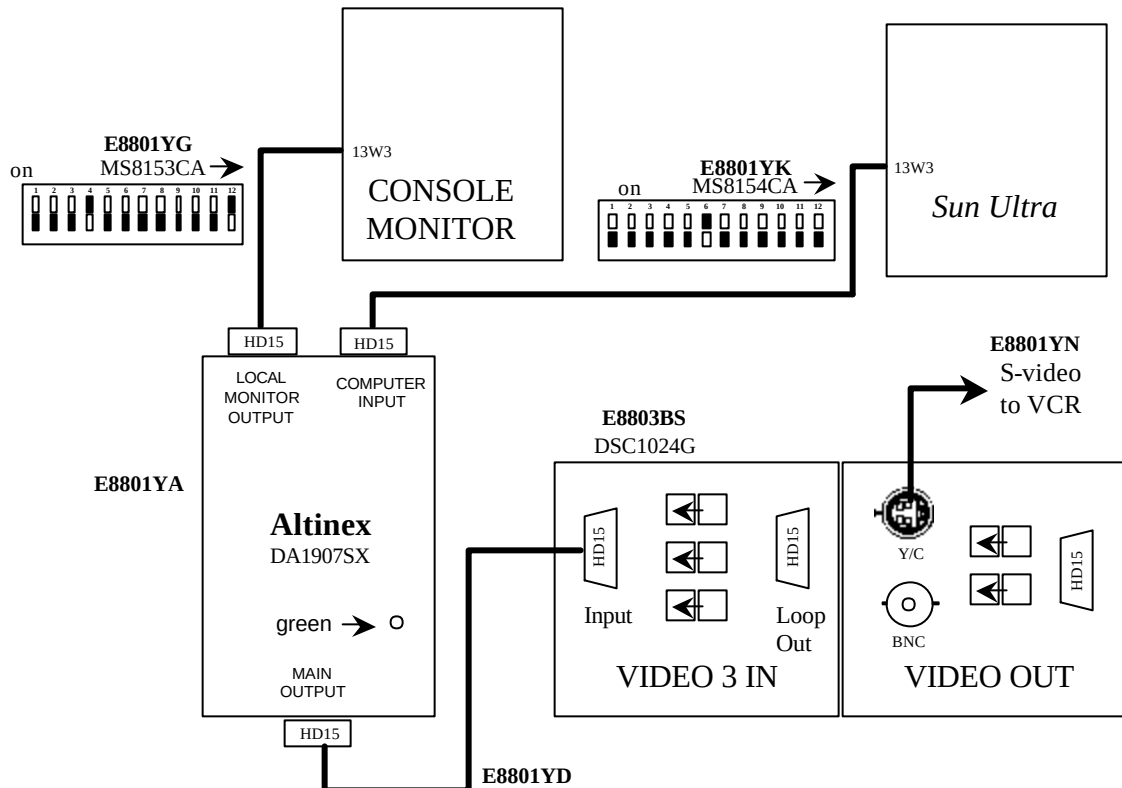


S9200LT MR Contour/Profile 3 (Color Video)

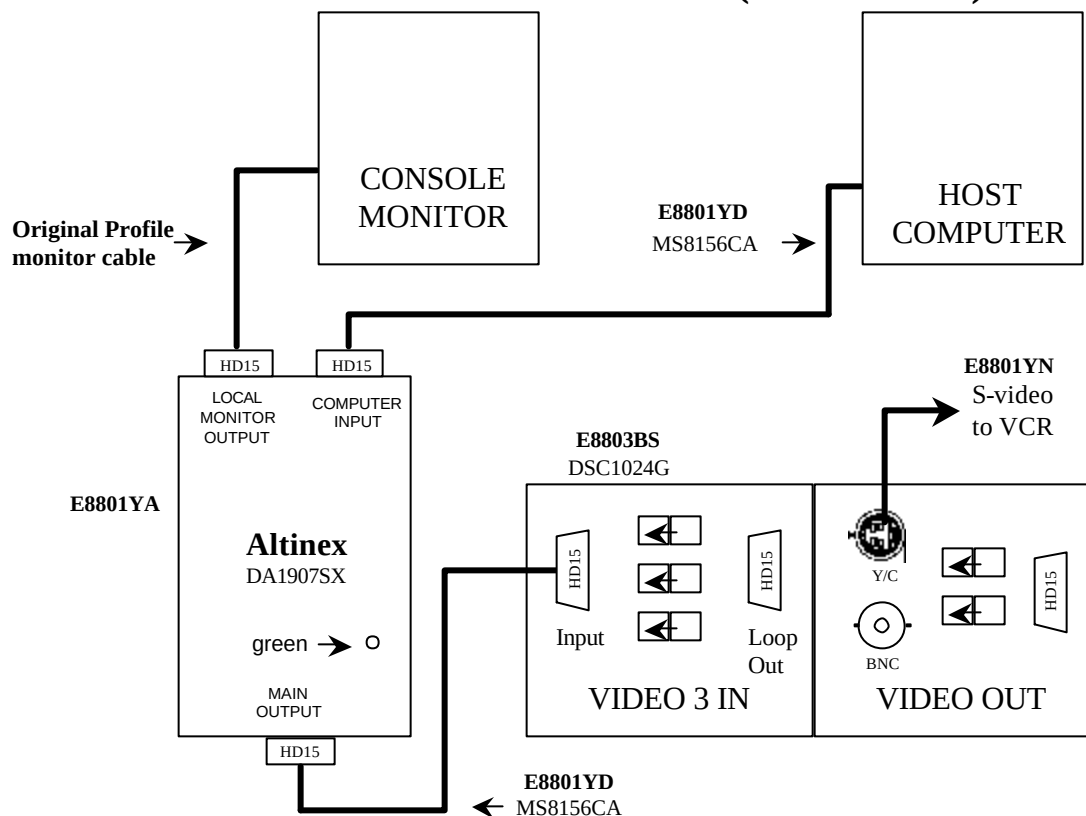

S9200LT MR Cardio/Neuro Vascular

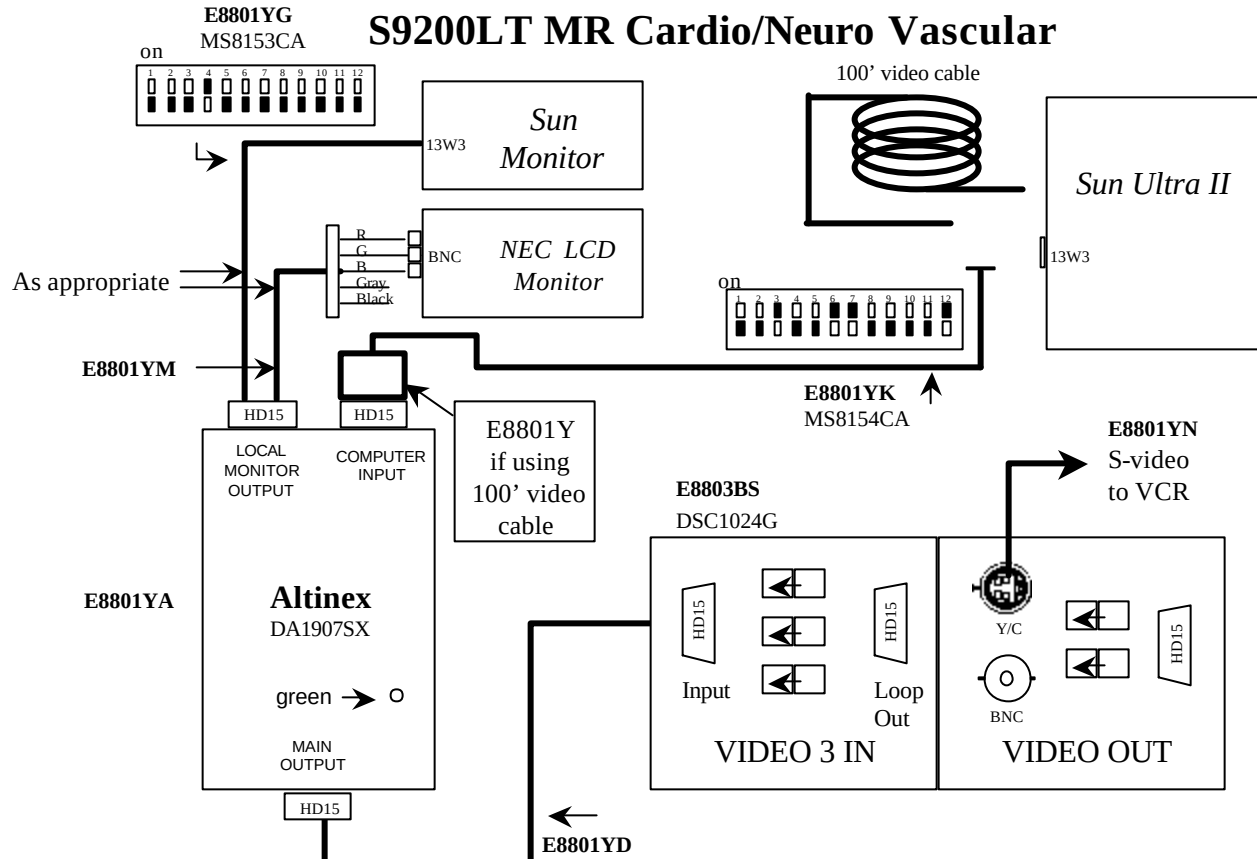


S9200LT MR Lightning LX

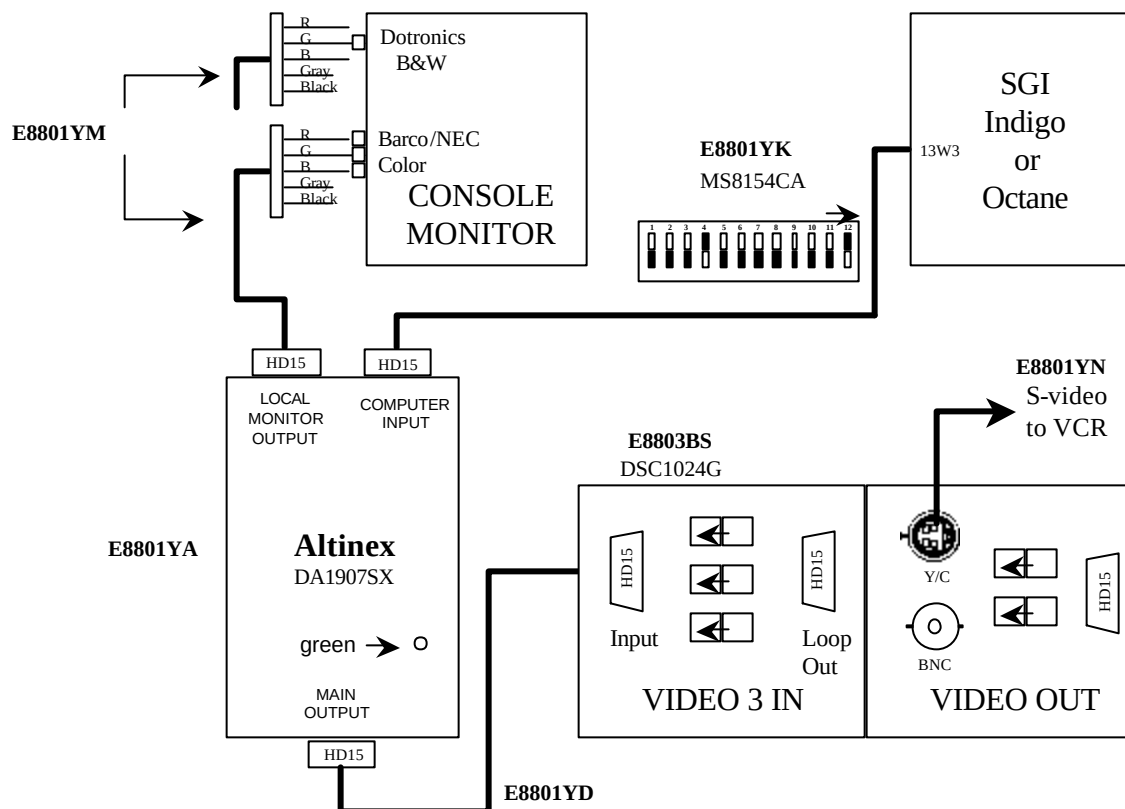

S9200LT Xray GEMNet

(with HP 725 Host Computer)

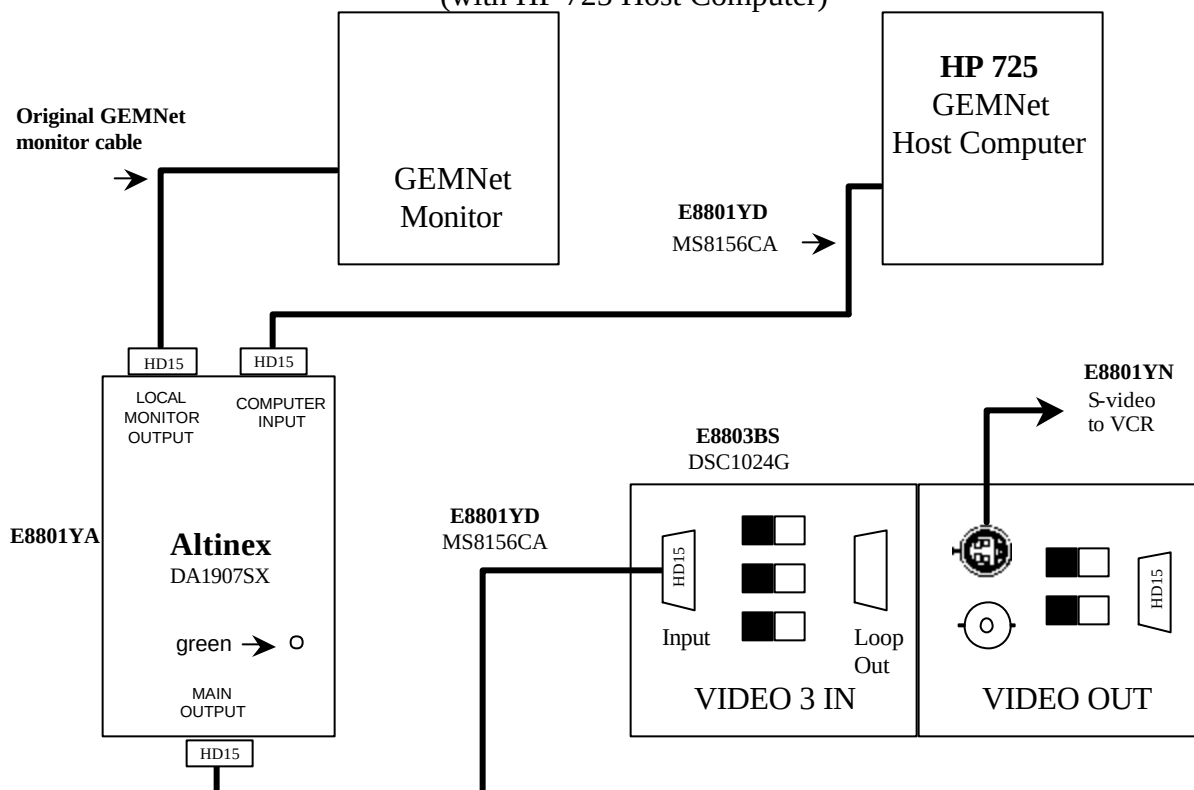


(with HP B180L Host Computer)




The diagram illustrates the E8801Y camera system configuration. It features a central **Altinex DA1907SX** unit. On the left, a box labeled **E8801YA** is connected to the Altinex unit via two **HD15** cables. The Altinex unit has a **LOCAL MONITOR OUTPUT** and a **COMPUTER INPUT** at the top, and a **MAIN OUTPUT** at the bottom. A **green** signal is shown entering the Altinex unit. To the right of the Altinex unit is a circle containing a **UP7300-MD Printer** with **Hi Scan Ports**. Below this circle is a box labeled **E7014PC UPZ-7000I**. A red note states: *Printer option not included in kit*. To the right of the printer is a box labeled **E8803BS DSC1024G** with **VIDEO 3 IN** and **VIDEO OUT** sections. The **VIDEO 3 IN** section has an **Input** and a **Loop Out**. The **VIDEO OUT** section has a **Y/C** output and a **BNC** output. A **HD15** cable connects the **VIDEO OUT** section to a box labeled **E8801YN** with **S-video to VCR** output.


S9200LT MR HP PC Linux Host

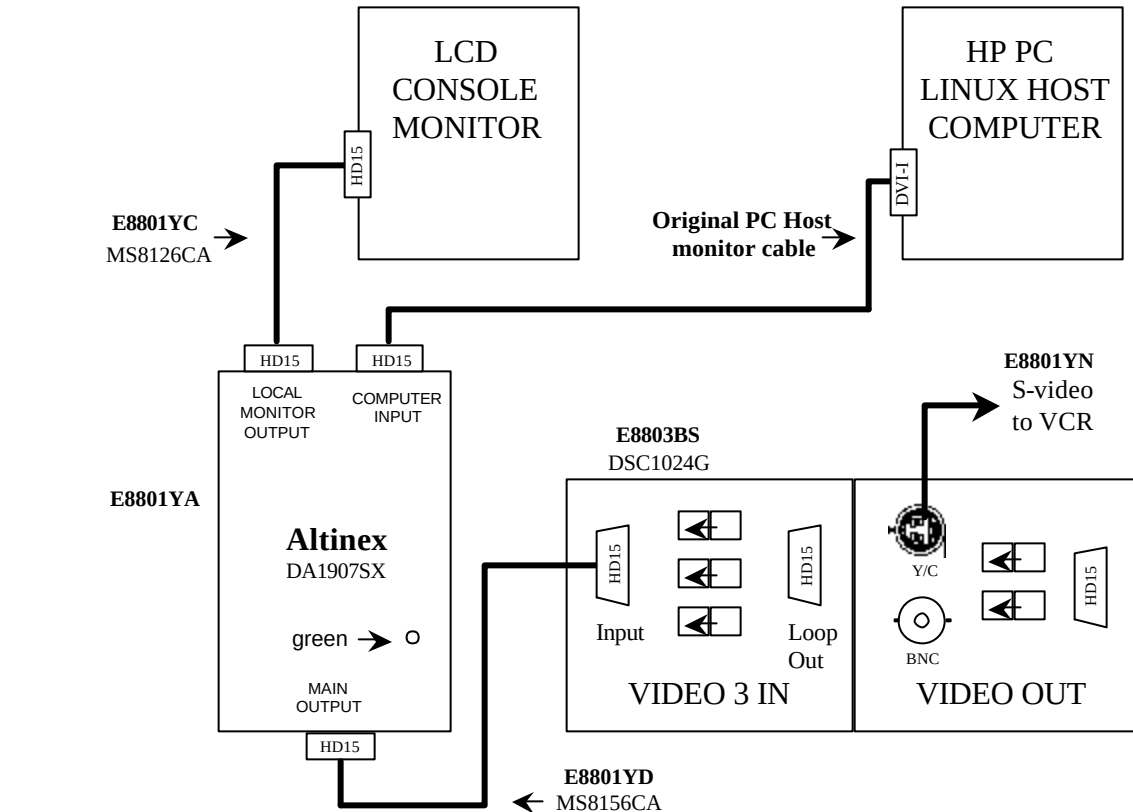

S9200LW MR Contour/Profile (Monochrome Video)

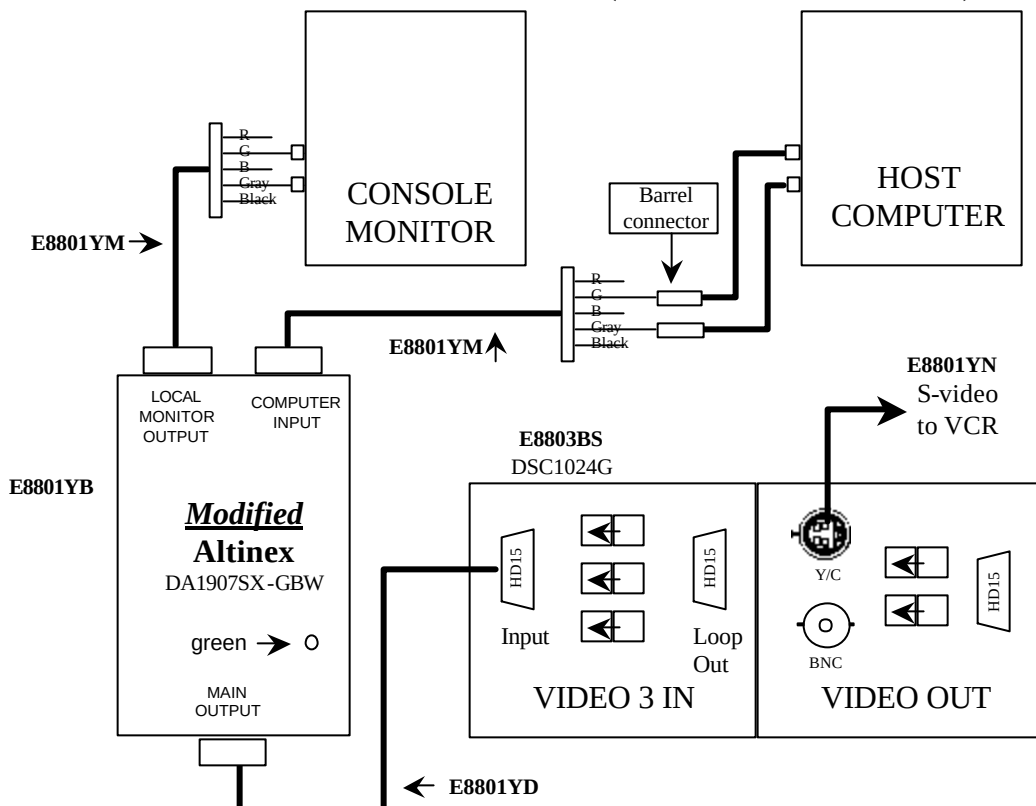

S9200LS Prestige/Advantx R&F, Vascular, Cardiac (DRS/DLX Monochrome Video)

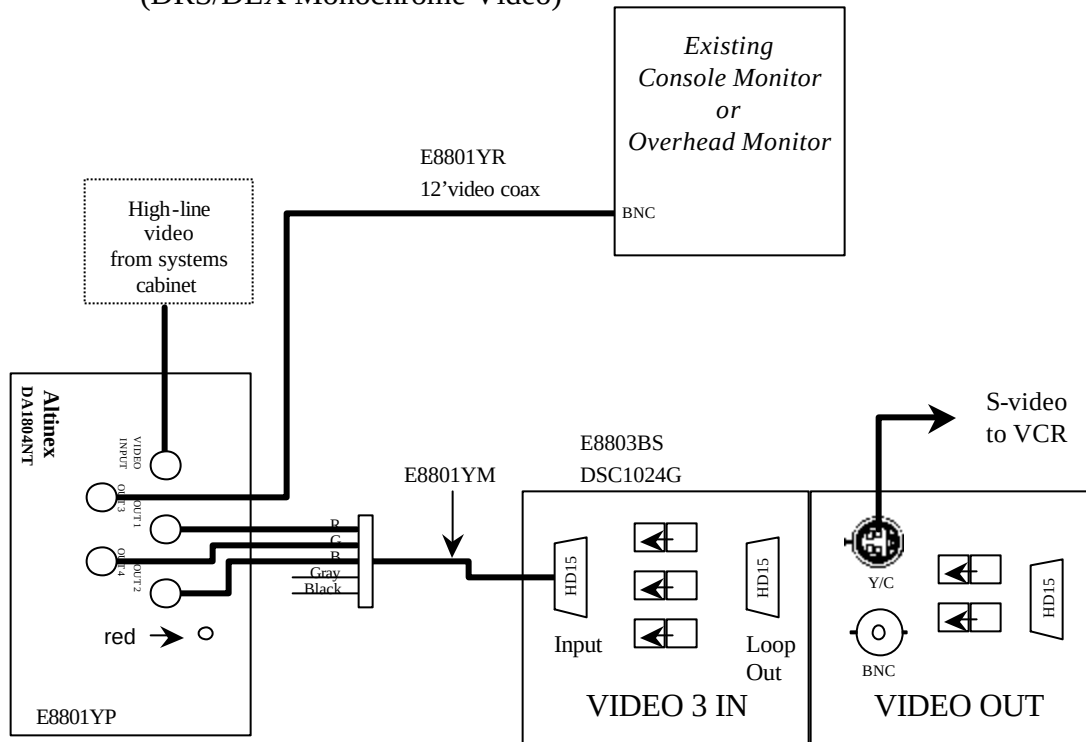

S9200LY MR Lightning LX

