GE Healthcare

Signa HDx Single Term Server Upgrade for RRI High Order Shim (HOS) Supply

Technical Publication

5177005
Rev 1

DAMAGE IN TRANSPORTATION

All packages should be closely examined at time of delivery. If damage is apparent, have notation “**damage in shipment**” written on **all** copies of the freight or express bill **before** delivery is accepted or “signed for” by a General Electric representative or a hospital receiving agent. Whether noted or concealed, damage **MUST** be reported to the carrier **immediately** upon discovery, or in any event, within **14** days after receipt, and the contents and containers held for inspection by the carrier. A transportation company will not pay a claim for damage if an inspection is not requested within this **14** day period.

File a report with

- Call 1–800–548–3366 and use option 8.
- Contact your local service coordinator for more information on this process.

Rev. 08/15/2003

© 2006, General Electric Company
GE Medical Systems, a General Electric Company, going to market as GE Healthcare
3200 N. Grandview Boulevard
Waukesha, Wisconsin 53188
USA

www.gehealthcare.com

Table of Contents

DESCRIPTION	PAGE
Table Of Contents	1
Language Policy For Service Documentation	3
Revision History	5
SECTION 1 – Overview	7
SECTION 2 – Upgrade Procedure	9

This page left blank intentionally.

WARNING

- THIS SERVICE MANUAL IS AVAILABLE IN ENGLISH ONLY.
- IF A CUSTOMER'S SERVICE PROVIDER REQUIRES A LANGUAGE OTHER THAN ENGLISH, IT IS THE CUSTOMER'S RESPONSIBILITY TO PROVIDE TRANSLATION SERVICES.
- DO NOT ATTEMPT TO SERVICE THE EQUIPMENT UNLESS THIS SERVICE MANUAL HAS BEEN CONSULTED AND IS UNDERSTOOD.
- FAILURE TO HEED THIS WARNING MAY RESULT IN INJURY TO THE SERVICE PROVIDER, OPERATOR OR PATIENT FROM ELECTRIC SHOCK, MECHANICAL OR OTHER HAZARDS.

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Revision History

Date	Rev	Change Description
5-04-2006	1	Initial Release

This page left blank intentionally.

SECTION 1 – OVERVIEW

1-1 OVERVIEW

This manual provides instructions to install a Single Term Server for 1.5T TwinSpeed or 3T94 systems equipped with a RRI type High Order Shim (HOS) supply. Plexus HOS supplies use CAN communication.

This Direction covers the following catalogs:

M3333SA – Single Term Server Kit Upgrade


This upgrade adds:

- Single Port Term Server
- Cables

1-2 SITE READY CHECK

Pre-installation work must be completed before equipment is delivered to avoid delays and confusion.

1-3 UPGRADE FLOWCHART

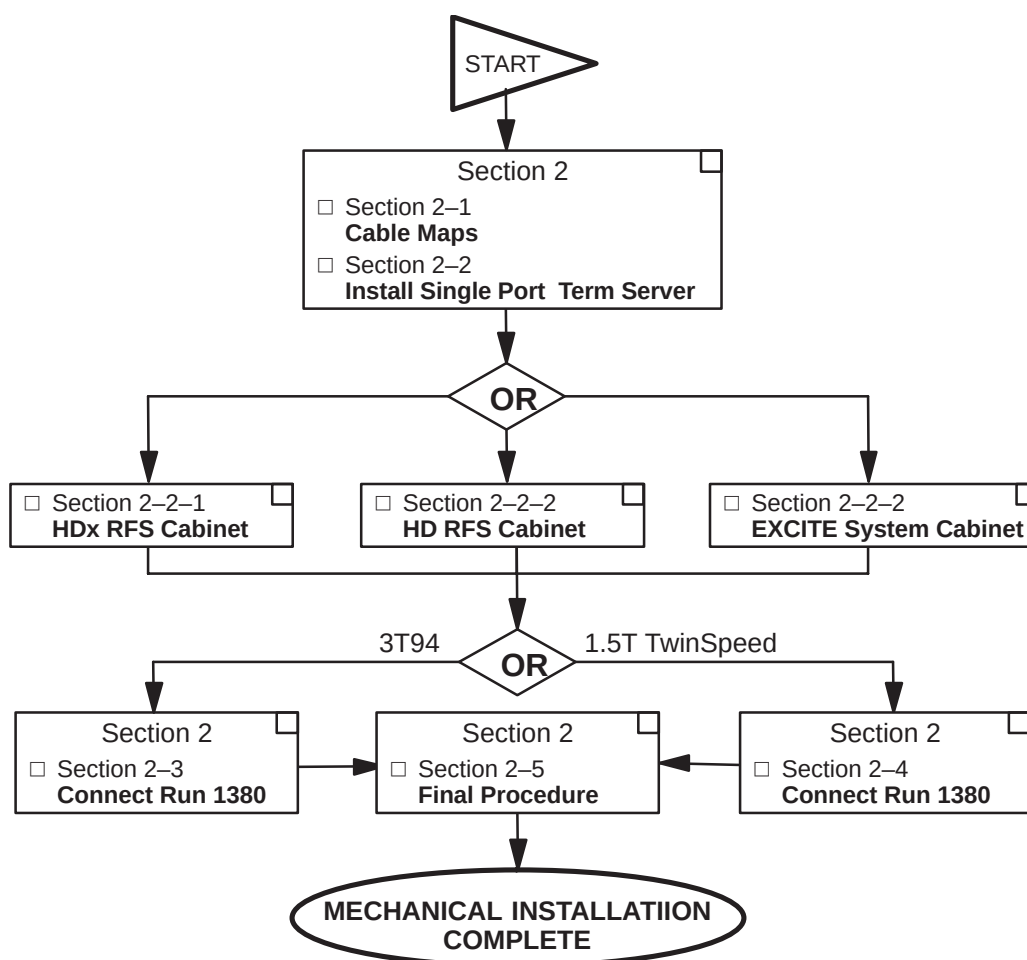1-4 FINAL COORDINATION

- . Replace all covers that have not been previously replaced.
- . Record and enter applicable data into applicable site configuration files and records.
- . Process product locator installation cards for all serialized components added to the system to include product locator card supplied with the upgrade serial number rating plate. Return to:

Product Locator File
P.O. Box 414, W-523
Milwaukee, WI 53201-0414

Note


Failure to fill out and return Product Locator Cards may result in failure of your site to receive future FMI's.

SECTION 2 – UPGRADE PROCEDURE


The following procedures are for installing the Single Term Server Upgrade Kit (5175780) in an HDx system. Refer to the applicable sections for a newly installed HDx RFS cabinet, or for an EXCITE system cabinet or HD RFS cabinet that was upgraded to HDx. There is also a section for installation of the upgrade in a 3T94 system with an Accessory cabinet.

2-1 CABLE MAPS

The cable maps below show the affected cables and components before and after upgrade for a 1.5T TwinSpeed system.

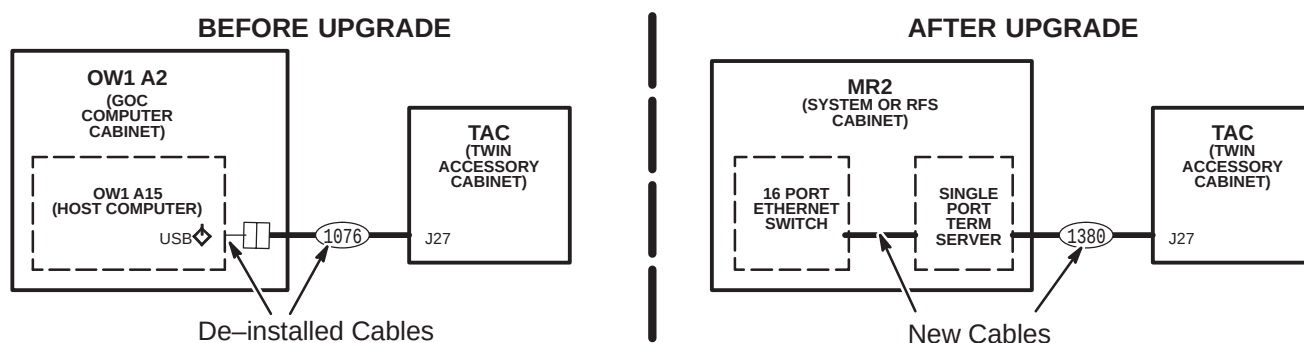The cable maps below show the affected cables and components before and after upgrade for a 3T94 system.

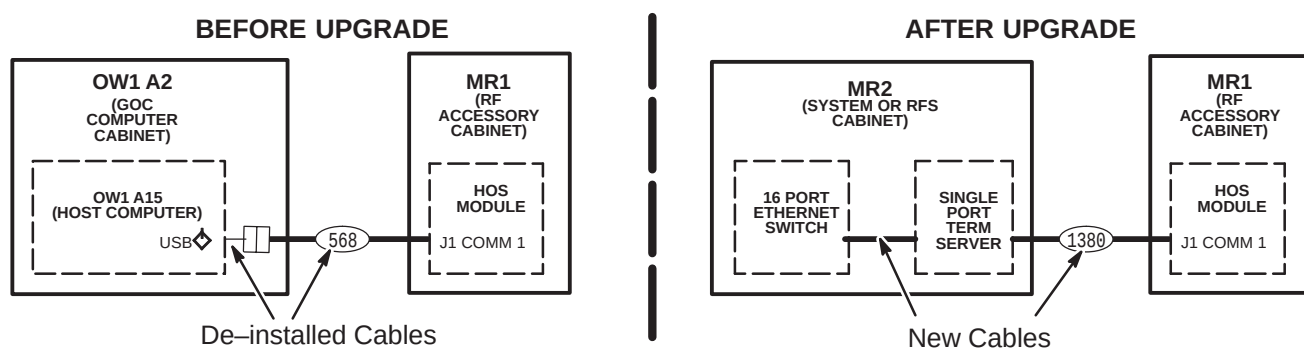2-2 INSTALL SINGLE PORT TERM SERVER

Of the following three sections, select the applicable cabinet for site being upgraded and perform installation of Single Port Term Server.

Note:


The Term Server may come in a box with an additional CD and cable (or terminator). Neither of these are used with this upgrade.

Before mounting the Single Port Term Server (5115759), record the 12 digit MAC address located on the side of the module. Write the address here:

_____ : _____


2-2-1 HDx RFS Cabinet

The following procedures are for installing this upgrade kit in an HDx RFS cabinet.

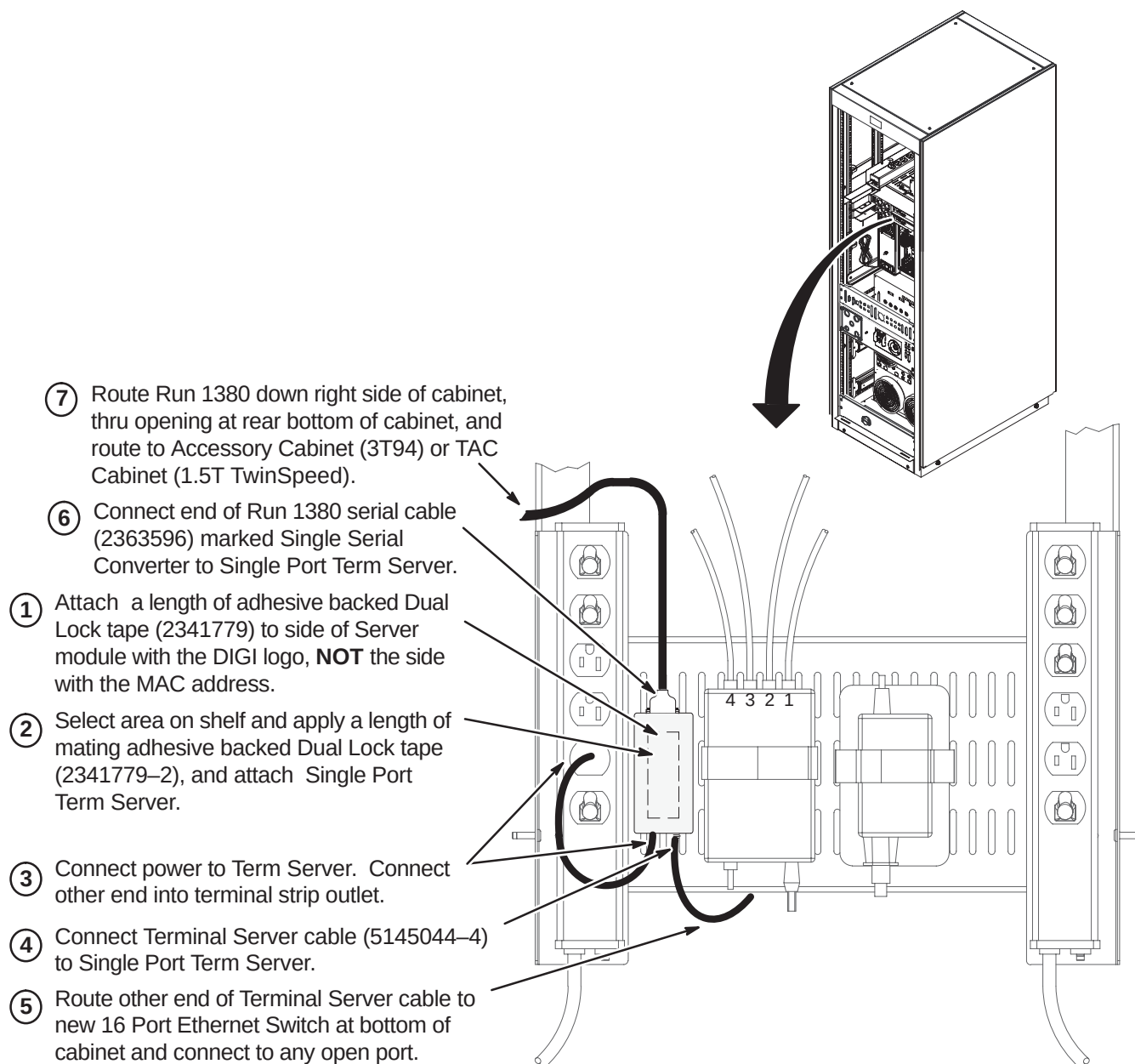
2-2-2 HD RFS Cabinet


The following procedures are for installing this upgrade kit in an HD RFS cabinet upgraded to HDx.

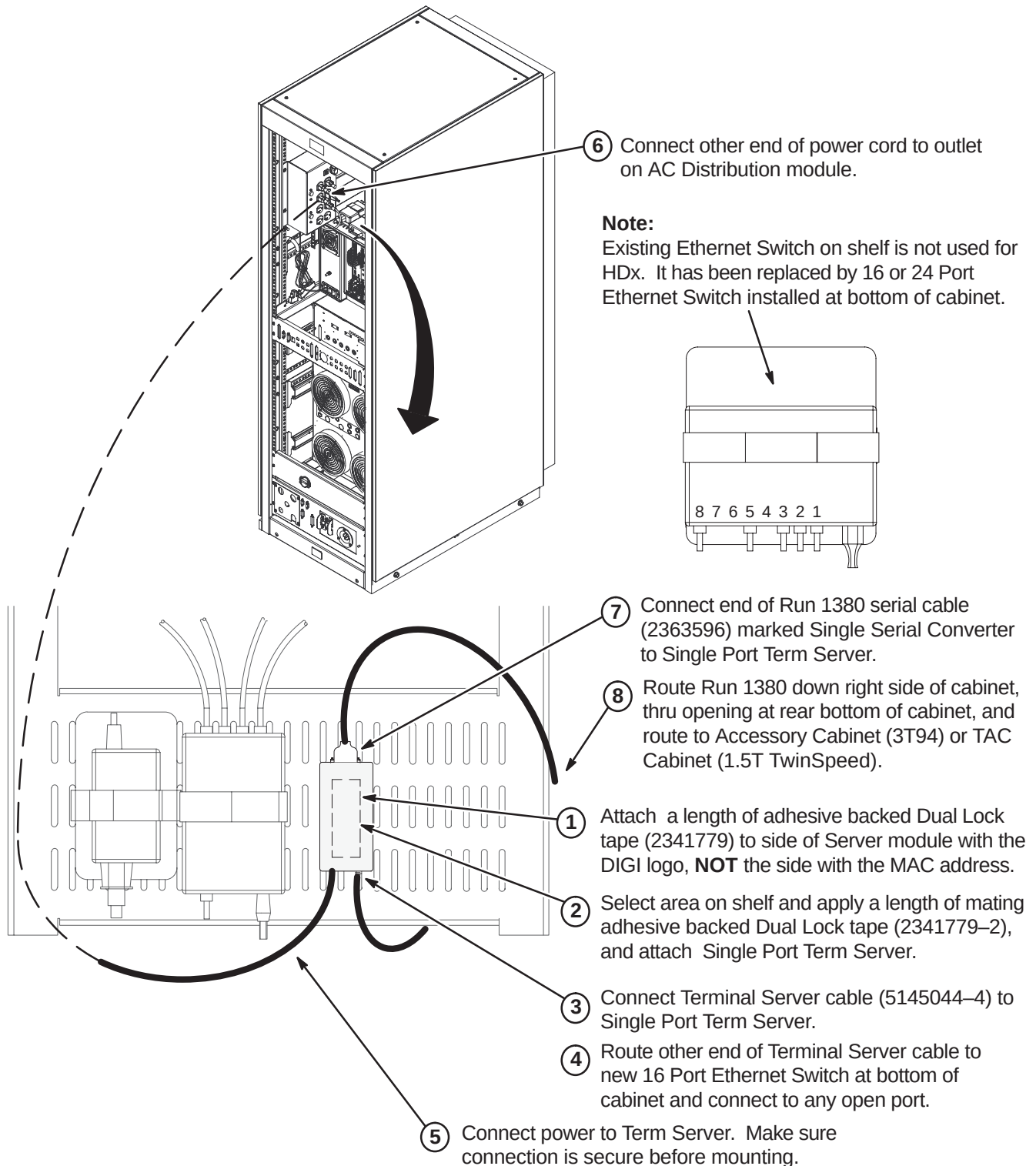2-2-3 SYSTEM CABINET

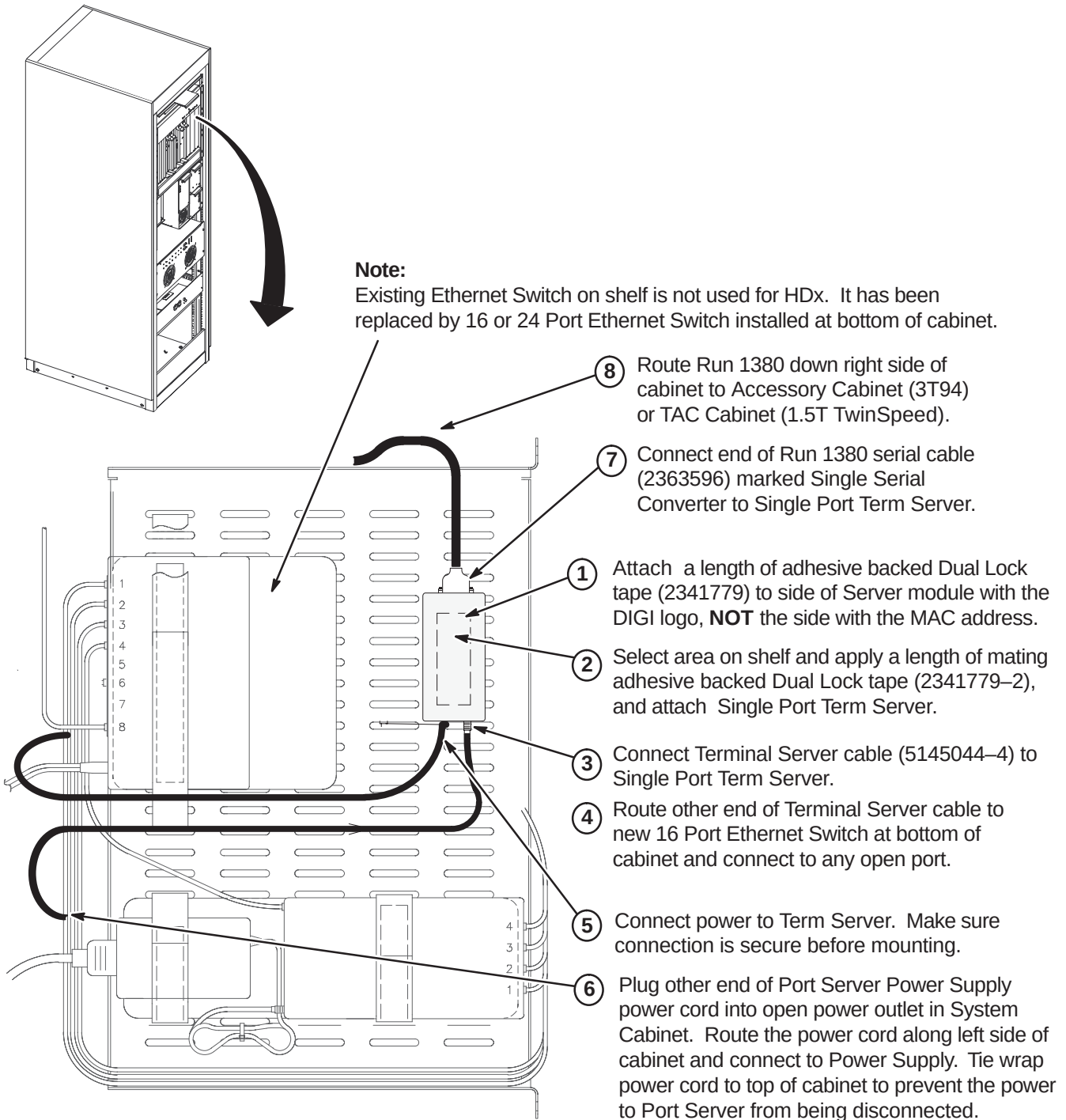2-3 CONNECT RUN 1380 TO HOS IN ACCESSORY CABINET

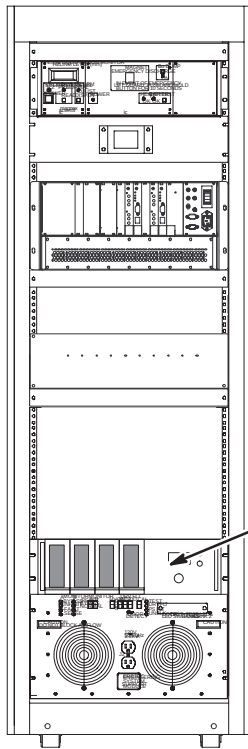**Note:**


Complete removal of obsolete cable is not required. It must, however, be disconnected at both ends as instructed.

- ① At rear of High Order Shim supply module in RF Accessory Cabinet, disconnect Run 568 from "J1 Comm 1" and remove cable.

Note: If there is anything else connected to J1 Comm 1, remove and discard.

- ② Other end of cable is connected to High Order Shim option cable in GOC. Disconnect and remove Run 568 at the GOC end.
- ③ Route and connect Run 1380 to J1 on HOS supply module.

2-4 CONNECT RUN 1380 TO TAC CABINET

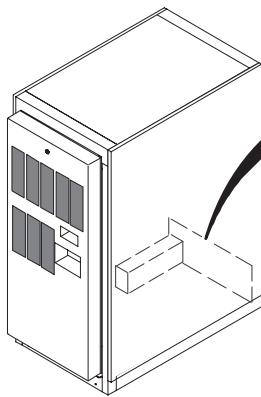- ① At rear of TAC Cabinet, disconnect Run 1076 from J27 and remove cable.

Note: If there is any other cable connected to J27, remove and discard.

- ② Other end of cable is connected to High Order Shim option cable in GOC. Disconnect and remove Run 1076 at the GOC end.
- ③ Route and connect Run 1380 from Single Term Server to J27.

2-5 FINAL PROCEDURE

1. Configure in Guided Install.
2. Test functionality.