



GE Medical Systems

Technical Publications

Direction 5133314

Revision 10

Signa[®] EXCITE HD 1.5T Field Replaceable Units (FRU)

Copyright © 2004 – 2009 by General Electric Company, Inc.

All Rights Reserved



GE Medical Systems

GE Medical Systems; Telex 3797371
P.O. Box 414, Milwaukee, Wisconsin 53201 U.S.A.
(Asia, Pacific, Latin America, North America)

GE Medical Systems – Europe; Telex 261794
Shortlands, Hammersmith, London W6 8BX U.K.



**COVER
BOTTOM
CROP LINE**

**COVER
BOTTOM
CROP LINE**

**5/16" RADIUS
CORNER TRIM
TYPICAL**



**COVER
LEFT
CROP LINE**

**COVER
FOLD
LINE**

**(2) CENTER STAPLES
FOR BINDING**

**COVER
RIGHT
CROP LINE**



**5/16" RADIUS
CORNER TRIM
TYPICAL**

SEND COMMENTS ABOUT THIS MANUAL TO:

GE MEDICAL SYSTEMS
MANAGER, SERVICE ENGINEERING, W-808
3200 N. GRANDVIEW BLVD.
WAUKESHA, WISCONSIN 53188
USA

This manual belongs to:

INTRODUCTION TO FIELD REPLACEABLE UNIT (FRU)

This direction is a guide to the FRU's.

Parts designated as a FRU are critical to maintain the safety and/or functionality of the system. FRU's are available on a Priority 1 basis from Service Parts Operation.

As a reference aid, the FRU's are listed under an upper-level part number where they are used. In most cases the upper-level part number is not a FRU. In some cases the parts listed under an upper-level assembly may have been nested in intermediate sub-assemblies which are not shown.

Where possible, FRU's which are part of Option Catalogs are shown under the area where they would be installed, rather than being listed in the Coils or Options Sections (SEE TABLE OF CONTENTS).

LIST OF EFFECTIVE PAGES

<u>PAGE</u>	<u>REV</u>	<u>PAGE</u>	<u>REV</u>	<u>PAGE</u>	<u>REV</u>
Cover Page .	10	11 to 14	7	28	9
1 to 2	10	15	8	29 to 40	7
3	7	16 to 17	7	41	8
4	8	18	9	42 to 44	7
5	6	19	10	45 to 46	8
6	5	20 to 21	7	47 to 48	7
7	7	22 to 23	8	49	8
8	9	24 to 26	7	50 to 51	7
9	7	27	8	52	10
10	10				

TITLE	PAGE
INTRODUCTION TO FIELD REPLACEABLE UNIT(FRU)	1
LIST OF EFFECTIVE PAGES	1
TABLE OF CONTENTS	2
PICTORIAL TABLE OF CONTENTS	3
MAGNET AND ENCLOSURE (ASSEMBLIES)	
BRM & TRM LIST (OF PARTS NOT SHOWN)	4
FRONT COVER, REAR END BELL	5
ENCLOSURE & VACUUM RELATED PARTS	6
FRONT & REAR ENDBELL, BRIDGE SUPPORT	6
SPLIT BRIDGE & MAGNET DOCKING PARTS	7
REAR PEDESTAL PARTS	8, 9
LOW PROFILE CARRIAGE ASSEMBLY	10
SRI/PAC & TOUCH-UP PAINTS	11
PATIENT SUPPORT TABLE	12, 13
CHILLER	15
PENETRATION PANEL	16
OPERATOR WORKSPACE	18
GOC COMPUTER CABINET	19
WORKSPACE SCIM AND KEYBOARDS	20
MAIN DISCONNECT PANEL	21
HFD-S/PDU CABINET	22
HFD/PDU CABINET	23
RFS CABINET	25
SYSTEM CABINET	26
MGD, RRF, & ASC MODULES	27, 28, 29, 30
SRFD CABINET	32
SRFD2 CABINET	33
TWIN ACCESSORY CABINET (TAC)	34, 35, 36
REPLACEMENT CABLES	39, 40
SERVICE TOOLS	41
COILS	43, 44, 45, 46
OPTIONS	48
SOFTWARE OPTION KEYS	49
BrainWave OPTIONS	50
REVISION HISTORY	52

INTRO TO FRU'S	PAGE 1
CHILLER PARTS	PAGE 15
CABLES	PAGE 39
COILS	PAGE 43
OPTIONS	PAGE 48
S/W OPTIONS	PAGE 49
BrainWave	PAGE 50

MAGNET ENCLOSURE:
SUBASSEMBLIES PGS
4, 5, 6, 8, 7, 10
TOUCH-UP PAINTS PAGE 11




2241999
 PATIENT AIR SUPPLY
 BLOWER CABINET (MG6)



PENETRATION
 PANEL
 PAGE 16


PATIENT
 SUPPORT
 TABLE
 PAGE 12

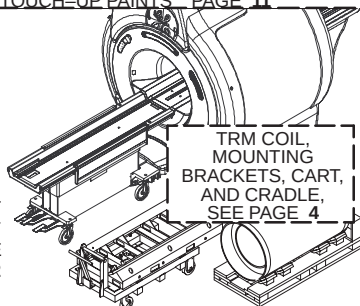

TRM COIL,
 MOUNTING
 BRACKETS, CART,
 AND CRADLE,
 SEE PAGE 4

MAGNET ROOM

OPERATOR/EQUIPMENT ROOM

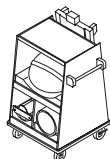

SERVICE TOOLS
 PAGE 41




2219400
 MAGNET
 MONITOR
 (MSM)




EXCITE OPERATOR
 WORKSPACE
 PAGES 18, 19, & 20




MAIN
 DISCONNECT
 PANEL
 PAGE 21

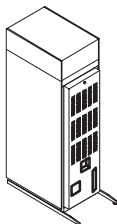


HFD-S/PDU
 PAGE 22
 HFD/PDU
 PAGE 23

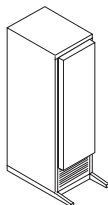

RFS CABINET
 PAGE 25
 SYSTEM CABINET
 PAGE 30

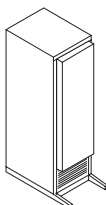

SRF CABINET
 PAGE 32
 SRFD2 CABINET
 PAGE 33

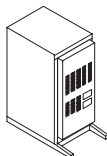

ACCESSORY
 CABINET
 PAGES 34
 THRU 36

PICTORIAL TABLE OF CONTENTS

NOTE: The following WideOpen & TwinSpeed common parts **Not Shown**:

5123812	BORE TEMPERATURE SENSOR, RUN 1260
5123813	OPERATOR DISPLAY CABLE TO SRI-3, RUN 1261
5133997	IIQ KIT, contains:
	BRIDGE SUPPORT FOR LCC MAGNET ENCLOSURES
	BRIDGE SUPPORT FOR S-SERIES MAGNET ENCLOSURES
	RADIAL SUPPORTS FOR ALL MAGNETS
5184777	RF SHIELD GROUND STRAP KIT
2324086-2	5 GALLON PAIL W/ETHYLENE GLYCOL WARNING LABEL

NOTE: The BRM parts **Not Shown**:

2379255	FOAM, GRAY, 116" LG. (Inserted between BRM and bore)
2297672	GLYCOL, 1 GALLON/3.8 LITERS
2124499-5	GRAD COIL LCC ONE-PIECE SHIELD
2144093-3	BRM GRADIENT COIL W/CART & CRADLE (Cx & S-series)
2144093-5	BRM-D GRADIENT COIL W/CART & CRADLE (LCC)
2241521	SPLICE KIT GRADIENT COIL SPLICE KIT
5110336	BRM GRADIENT CABLE LEAD BOARD KIT
5127441-2	GRADIENT CABLE CLAMP BLOCK KIT
2208999-3	SHORT RF COIL & CASE, Order for RF COIL Replacement (LCC)
2208999-6	LONG RF COIL & CASE, Order for RF COIL Replacement (Cx & S-series)
5113879*	RF SHIELD FOR BRM COIL

NOTE: The following TwinSpeed 1.5T FIXED W/TRM & Enclosure related parts **Not Shown**;
MISCELLANEOUS:



2297760	VACUUM GREASE, MSDS# 8366442
2203803	SPLIT O-RING
2242756	RELIEF VALVE, 3/8 NPT X 7 PSIG, GEN
2299843	MAGNET PRESSURE TRANSDUCER, SUNX DP2-41N
46-294231G1	MAGNET RUNDOWN UNIT
46-294231P1	BUSSMAN, F1, 110VAC .400 GDC
46-264758G1	REMOTE INTERCOM BOARD, MG3 A2 A1

TRM:


2297672	GLYCOL, 1 GALLON/3.8 LITERS
5113878*	RF SHIELD FOR TRM COIL
2144093-7	TRM WITH MOUNTING BRACKETS, CART, AND CRADLE
2218148	COLD SHRINK (FOR REAR ENDBELL AND TRM COIL)
2208999-8	RF BODY COIL AND CASE
2299498	COOLING CIRCUIT REPAIR KIT
2303520	DIODE SLOT BOARD
5127441	GRADIENT CABLE CLAMP BLOCK KIT
2293612	MALE CONNECTOR 3/4" HOSE TO 3/4" NPT

***2389049** RF Shield Change Out Kit – MUST be ordered for the RF Shield replacement

4 MAGNET AND ENCLOSURE (ASSEMBLIES)

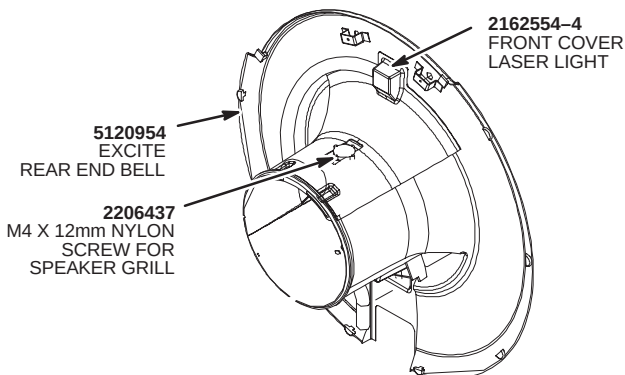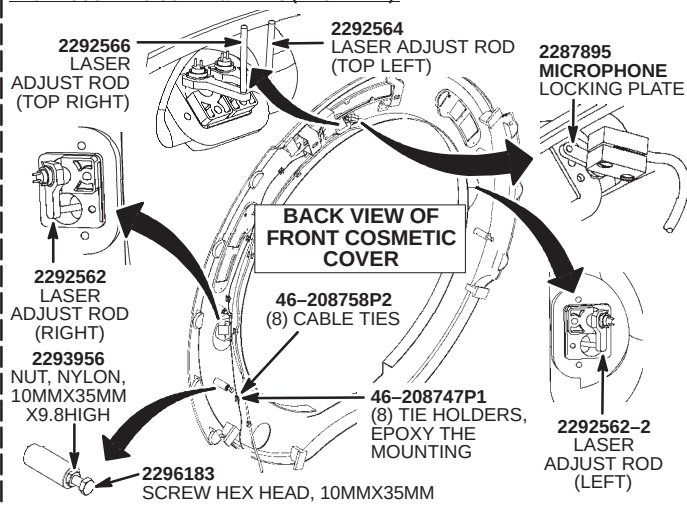
FRONT COSMETIC COVER & PARTS (BACK VIEW):

ENCLOSURE & VACUUM RELATED PARTS FOR TRM:

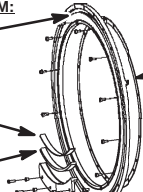- 2293323**
 O-RING BETWEEN ENDBELL
 AND ALUMINUM RING
2293771
 RADIAL SHIM SPACER, 1.5MM
2293771-2
 RADIAL SHIM SPACER, 0.8MM
2241176
 RUBBER ISOLATION PAD
- 2246966**
 MAGNET
 RING
 GASKET
- 

**FRONT ENDBELL PARTS:
(TRM ONLY)**


- 2295639**
 GASKET,
 ENDBELL
 O-RING, TO RF TUBE
 ACCESS **2294358**, 4MM
 PANELS **2294358-2**, 6MM

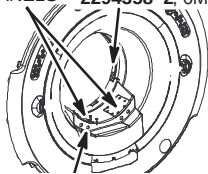

- 2297761**
 O-RING, WATERLINE
 FEEDTHROUGH

**COIL ASSEMBLY:
(TRM ONLY)**

- 2131953-5**
 MOUNTING PAD,
 1/4 INCH THICK
 NEOPRENE

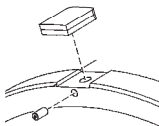

- 2293535**
 SET SCREW,
 MOUNTING PAD


**BRIDGE SUPPORT PARTS:
(TRM ONLY)**

- 2293539**
 BRIDGE
 SUPPORT
 BUSHING

- 2293956**
 NYLON HEX NUT
 M10 X 1.5MM
 X 9.5MM HIGH

**REAR ENDBELL ASSEMBLY:
(TRM ONLY)**

- 2284170**
 REAR ENDBELL
 VACUUM ASSEMBLY;
 & **2218148**
 COLD SHRINK
 (FOR REAR ENDBELL &
 TRM COIL)
 & **2218149**
 COMPRESSION
 SPLICER
 (FOR REAR ENDBELL)

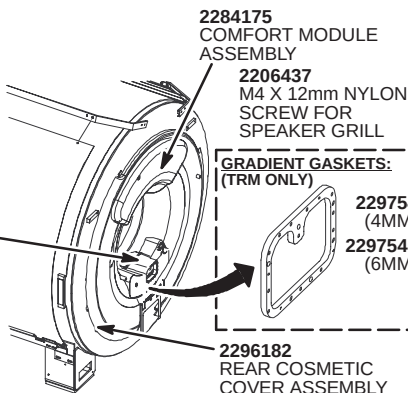

- 2284175**
 COMFORT MODULE
 ASSEMBLY


- 2206437**
 M4 X 12mm NYLON
 SCREW FOR
 SPEAKER GRILL

**GRADIENT GASKETS:
(TRM ONLY)**

- 2297540**
 (4MM)


- 2297540-2**
 (6MM)

- 2296182**
 REAR COSMETIC
 COVER ASSEMBLY

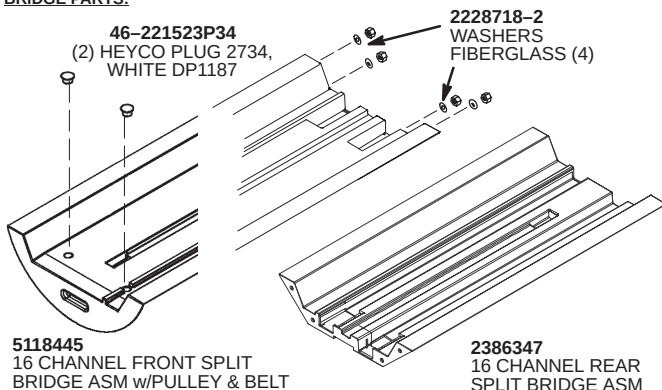BRIDGE PARTS:


NOTE: The following part is **Not Shown**.



2386348	FRONT BRIDGE FIBERGLASS
5117304	FRONT BRIDGE PULLEY ASM
2293538	BRIDGE EDGE FOAM
46-282436P1	CARRIAGE TRAVEL SENSOR FLAG

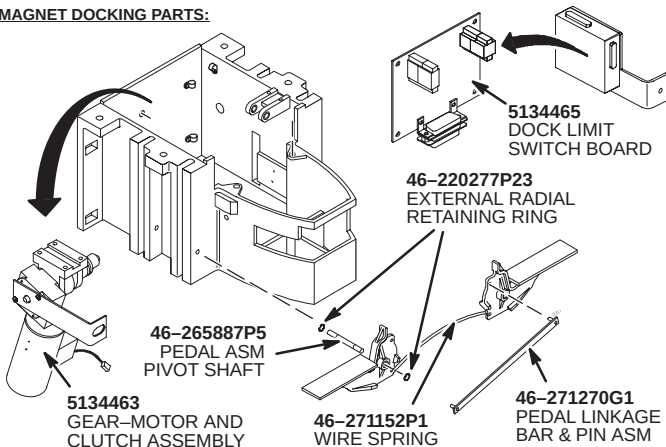
MAGNET DOCKING PARTS:**MAGNET AND ENCLOSURE (ASSEMBLIES)**

REAR PEDESTAL

Note: See next page for Hybrid Splitter information

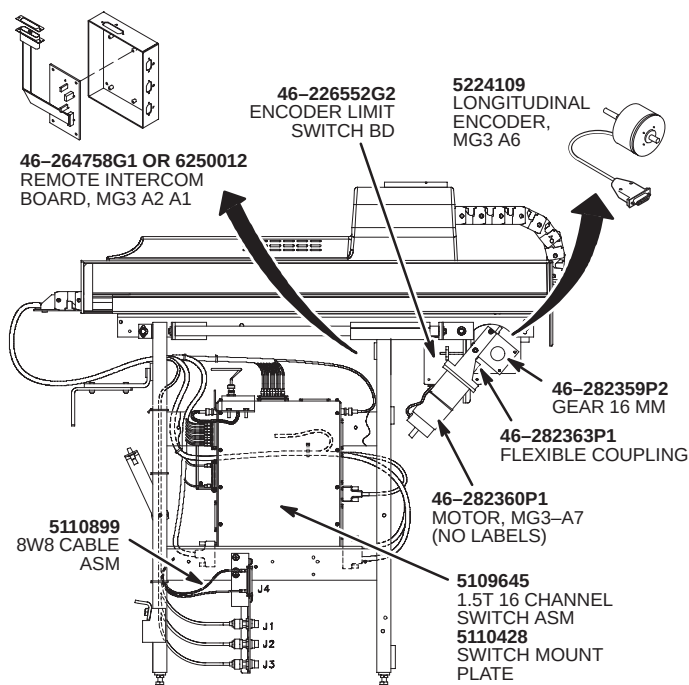


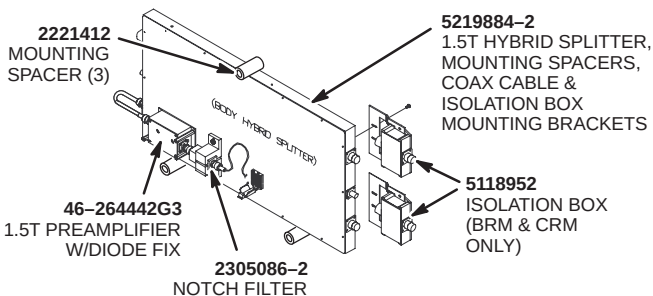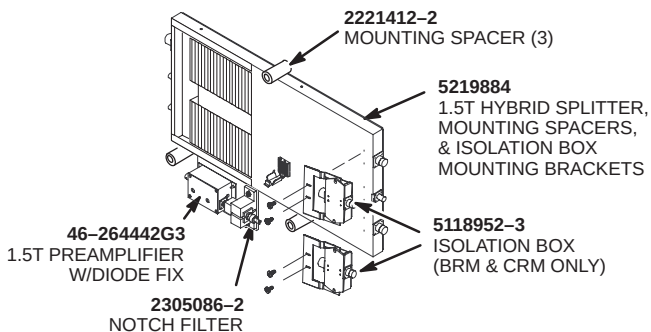REAR PEDESTAL**1.5T HYBRID SPLITTER – TELEDYNE:****1.5T HYBRID SPLITTER – NETCOM:**

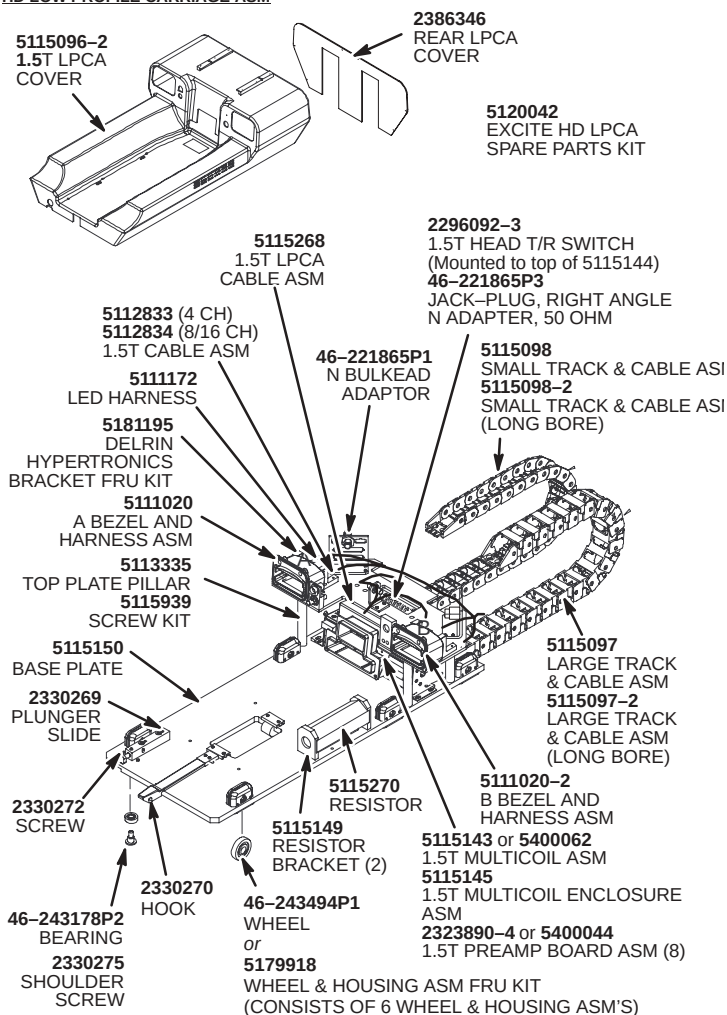
HD LOW PROFILE CARRIAGE ASM

SRI/PAC PLATFORM:

5112668
NON-LCC
PAC-2B MODULE

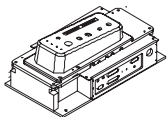

5122500
ECG LEAD KIT

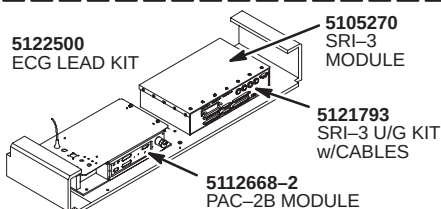

5105270
SRI-3
MODULE

5121793
SRI-3 U/G KIT
w/CABLES


5112668-2
PAC-2B MODULE

TOUCH-UP PAINTS:

2183491-3 GRAY #1 0.6 OZ. NON-AEROSOL

46-303460P3 STEEL BLUE 0.6 OZ. NON-AEROSOL


46-303460P4 DARK STEEL BLUE 0.6 OZ. NON-AEROSOL

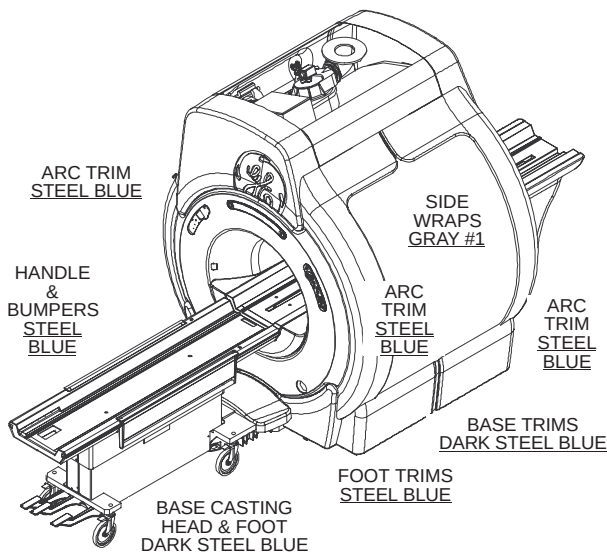


NOTE: To obtain enclosure part numbers not found in this manual, refer to Direction 2138247, *Signa Illustrated Parts List*, for enclosure covers, endbells, and other parts.

PATIENT TABLE – INFINITY and HDx

NOTE: Refer to page 13 for information on bumper FRU's for HD, Infinity, and Pre-Infinity patient tables.

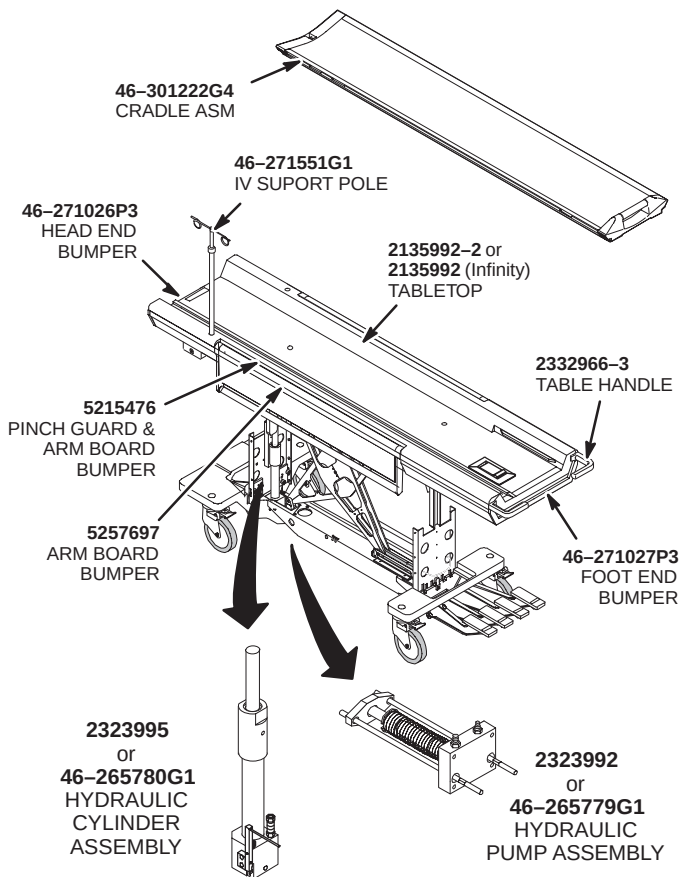

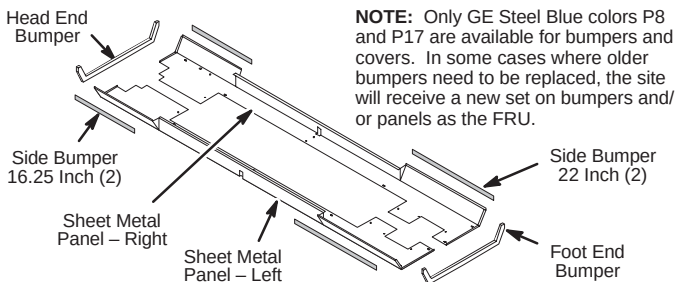PATIENT TABLE – BUMPERS AND PANELS


INFINITY EXISTING PART	INFINITY REPLACEMENT FRU
46-271022G1, 2393106, Panel – Left 46-271021G1, 2393106-2, Panel – Right	5179706 , New panels and side bumpers FRU kit for Infinity table
5154908-4, Panel – Left 5154908-3, Panel – Right	5179706 , New panels and side bumpers FRU kit for Infinity table
46-271025P5, 22 Side Bumper	5161172-3 , Side Bumper 22, Color P17
46-271025P6, 16.25 Side Bumper	5161172-4 , Side Bumper 22, Color P17
46-271027P3, Foot End Bumper	46-271027P3 , Foot End Bumper
46-271026P3, Head End Bumper	46-271026P3 , Head End Bumper
Pre-INFINITY EXISTING PART	Pre-INFINITY REPLACEMENT FRU
5154908-4, Panel – Left (Grey) 5154908-3, Panel – Right (Grey)	5179706 , New panels and side bumpers FRU kit for Infinity table
22 Side Bumper (Accent Grey or Blue)	5161172-3 , Side Bumper 22, Color P8
16.25 Side Bumper (Accent Grey or Blue)	5161172-4 , Side Bumper 16.25, Color P8
Foot End Bumper (Accent Grey or Blue)	46-271027P3 , Foot End Bumper
Head End Bumper (Accent Grey or Blue)	46-271026P3 , Head End Bumper

[illegible]

SYSTEM COOLING CABINETS:

REFER TO ELLIS & WATTS SERVICE
MANUALS FOR FRU INFORMATION

REFER TO AIRSYS
SERVICE MANUALS
FOR FRU INFORMATION

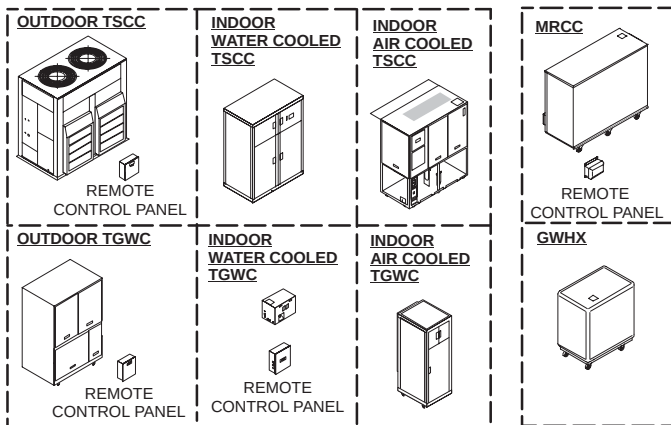

2343325 AIRSYS HEAT EXCHANGER (GWHX)

2341077-2 AIRSYS QUICK DISCONNECT KIT

2341077-3 AIRSYS REMOTE CONTROL PANEL

2230235 E & W MODEL LTL-4 MOBILE HEAT EXCHANGER

2230235-3 E & W PUMP/MOTOR SET

2230235-12 E & W LEVEL SWITCH

2230235-24 QUICK-DISCONNECT FACILITY

2230235-25 QUICK-DISCONNECT COIL

2298100-4 E & W MODEL LTL10 HEAT EXCHANGER CABINET

2298100-5 CONTROL PANEL FOR LTL10 ONLY

2297672 GLYCOL USED WITH E & W OR AIRSYS

NOTE: The following parts are **Not Shown**.

2222564-22 LYTRON CHILLER, 50Hz OR 60Hz

2222564-23 LYTRON CIRCUIT BREAKER

2222564-24 LYTRON PUMP

2222564-25 LYTRON MOTOR

2222564-26 LYTRON KEY PAD


2222564-27 LYTRON FLOW SWITCH

2222564-28 LYTRON HOSE FITTINGS


2222564-29 LYTRON PUMP with NYLON COUPLING

2222564-31 LYTRON NYLON MOTOR COUPLING

2260427 COOLANT (1 GALLON) (LYTRON ONLY)

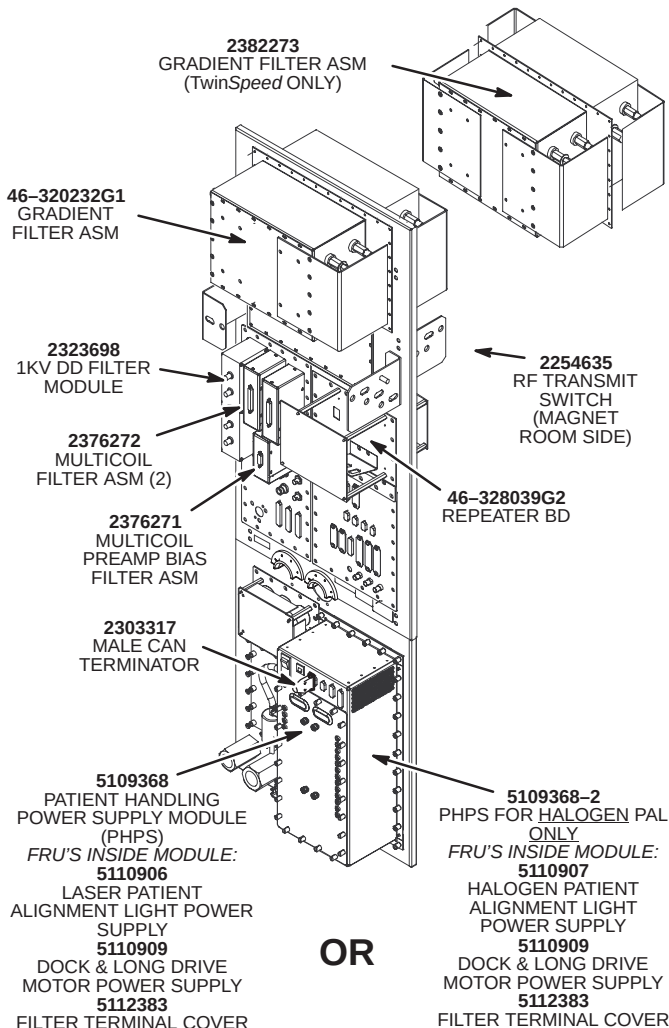

[illegible]

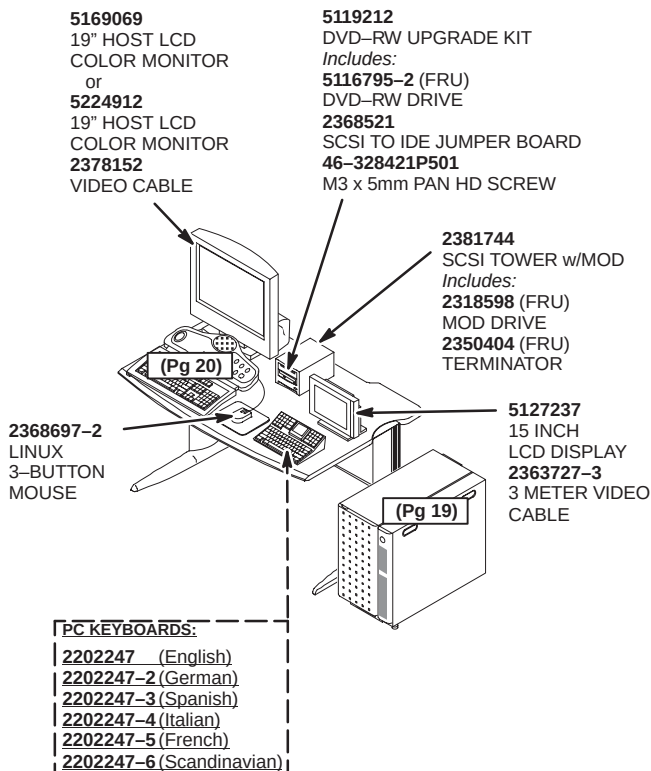

NOTE: The following parts are Not Shown.


2156158-50 S/W COLLECTOR FOR EXCITE HD (12.x) M4

2389540 SCAN TOOLS MOD


2376768 GAMMA PATCH CD-ROM FOR LCD MONITORS

GOC COMPUTER CABINET WITH HP8200:**NOTES:**


- The HP8200 has visible USB ports located on the right side of front of PC cover.
- The HP8200 cannot be replaced with the HP8000.
- If substituting HP8200 for HP8000, must order 2378152 video cable for waveform LCD monitor

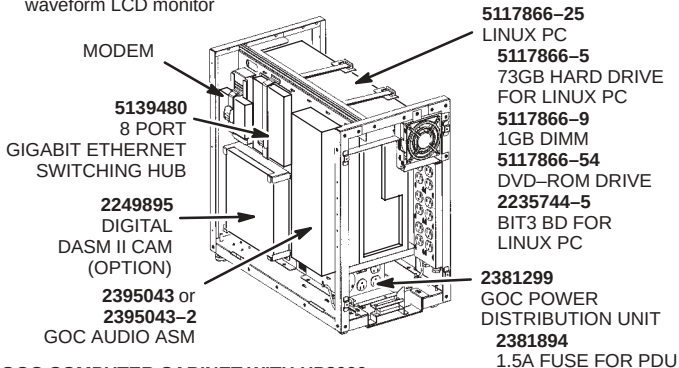**GOC COMPUTER CABINET WITH HP8000:****NOTE:**

- The HP8000 can be replaced with the HP8200.

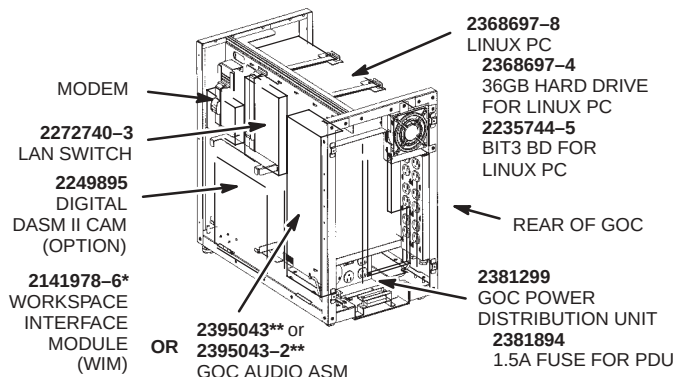**NOTES:**

- * If WIM is not available, order part 2141978-5 and replace cable between motherboard J1 and CIF J3.
- ** If replacing WIM with GOC AA, order cable 2395735.

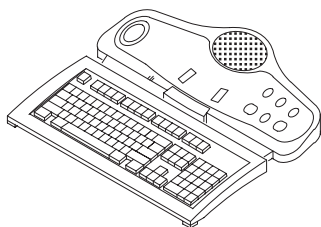LATEST VERSIONS:**KEYBOARDS ONLY:**

2275756 (English)
*2275756-2 (French)
*2275756-3 (German)
*2275756-4 (Scandinavian)
*2275756-5 (Italian)
*2275756-6 (Portuguese)
*2275756-7 (Spanish)

*** NOTE:**

These types of Keyboard/SCIM Assemblies were not available at the time this book was released, see below for earlier types.

NOTE: The following parts are Not Shown.

2262158-9	CABLE, 50 PIN SCIM KEYBOARD
2289697-2	SCIM ^W /ENGLISH OVERLAY
2289697-8	MR SCIM GLOBAL OVERLAY SET

NOTE: Refer to the Directions on CD-2160623 (Signa Service Methods) for additional product information.

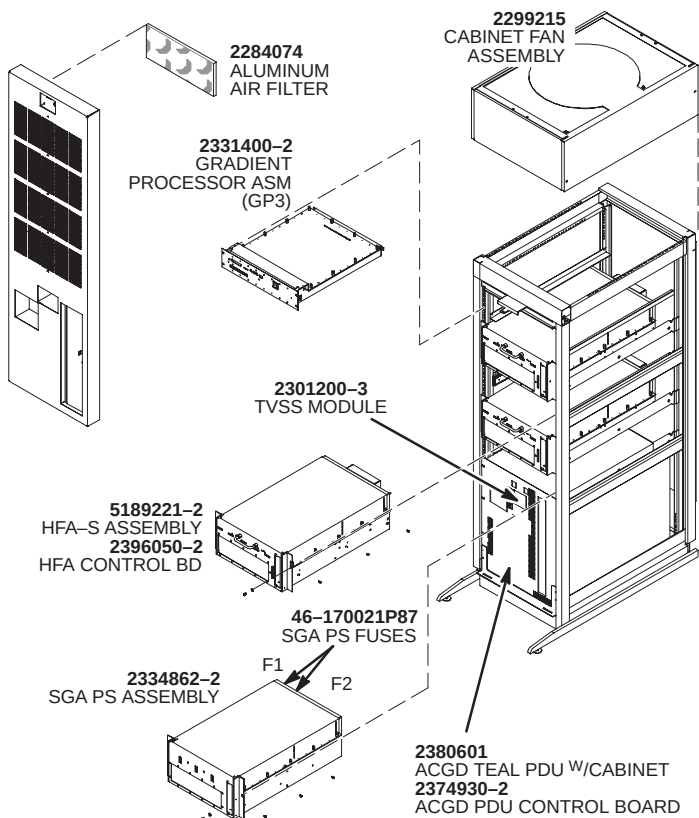


TABLE 1-12

LOCATION	AMPS	VOLTS	GEMS PART NUMBER	MANUFACTURER	
				BUSS	G-S
FUSE F1	1.25	480	2276547-2	SNQ-R-1-1/4	ATQR1-1/4
FUSE F2	3	120	2276547-3	FNQ-R-3	ATQR3
FUSE F3	1	12VDC	2276547-8	LP-CC-1	ATDR1
FUSE F4	0.12	480	2276547-9	SNQ-R-8	ATQR1/8
FUSE F5	125	600	2276547-19	LPJ125SP	AJT125
FUSE F6	50	600	2276547-20	LPJ-50SP	AJT50
FUSE F7	20	600	2276547-22	LPJ-20SP	AJT20
FUSE F8	5	600	2276547-21	FRS-R-5	TRS5R
FUSE F9	8	400	2276547-6	FRS-R-8	TRS8R
FUSE F10	10	120	2276547-11	FNQ-R-10	ATQR10
FUSE F11	10	120	2276547-11	FNQ-R-10	ATQR10


2276547-14 12 VDC BATTERY

2276547-16 12 VDC SUPPLY

HFD-S/PDU CABINET ASSEMBLY W/SGA/SGAPS/GP


NOTE: The following Part is Not Shown.

2303317 CAN Terminator

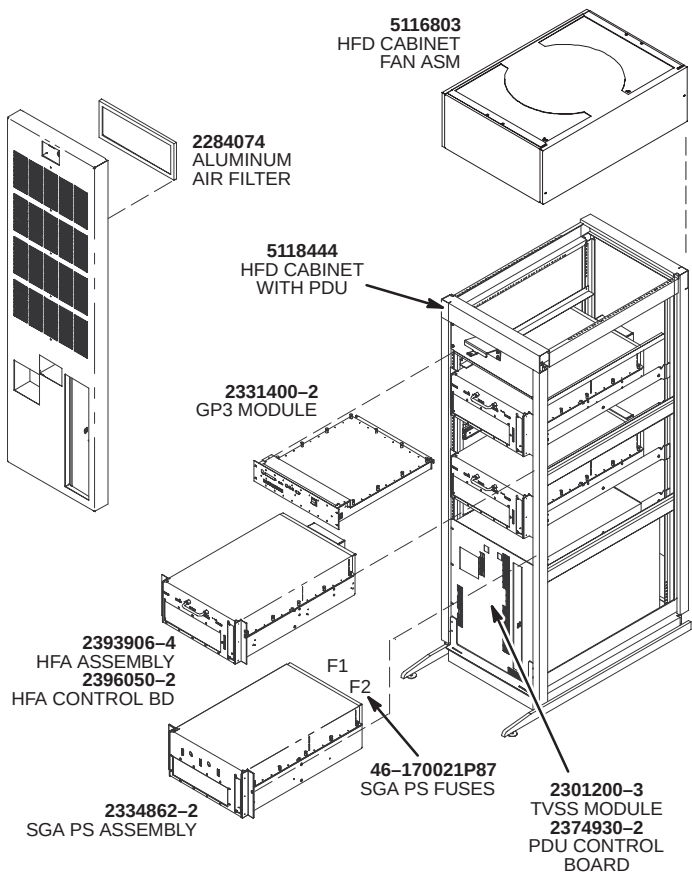HFD/PDU CABINET

24


NOTES

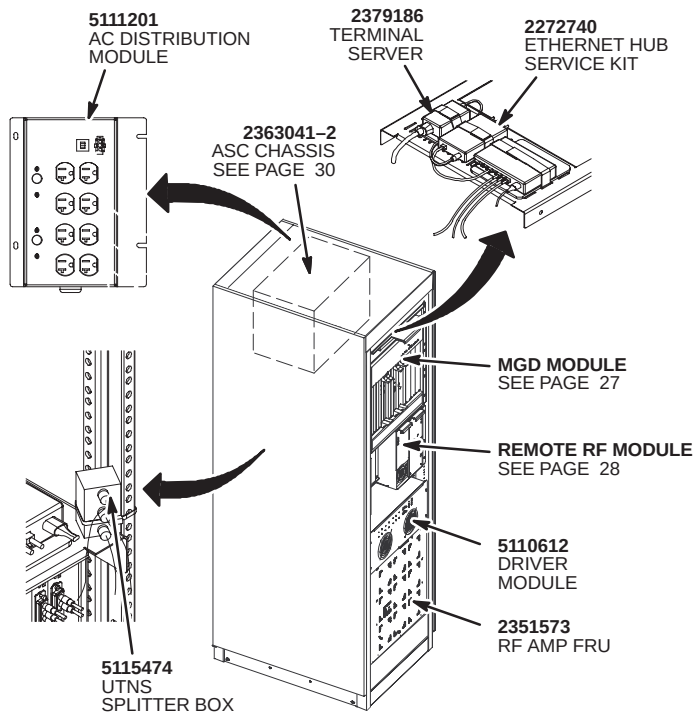

NOTE: The following items are Not Shown.

5117087 RRF SERVICE DIAGNOSTIC KIT


46-296817P1 50 OHM BNC TERMINATOR

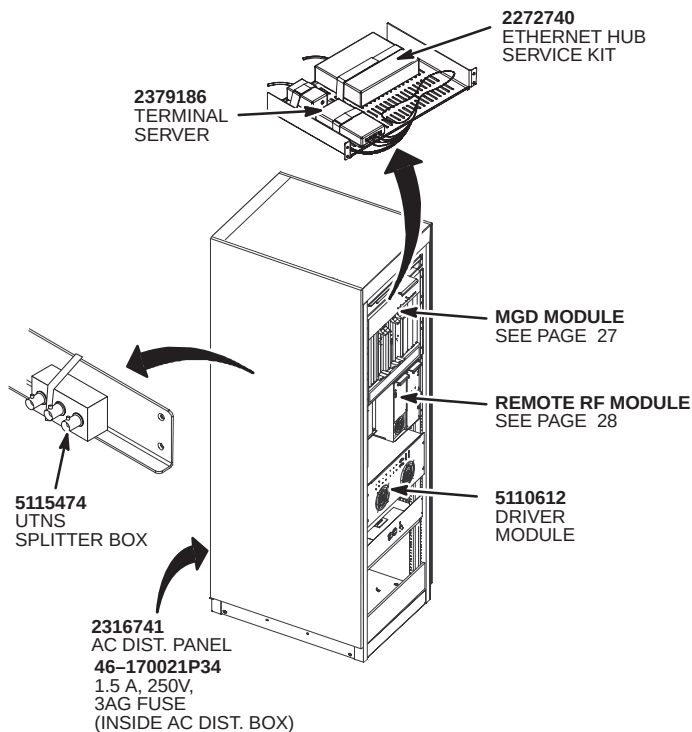EXISTING SYSTEM CABINET UPGRADED TO EXCITE HD

NOTE: The following items are Not Shown.

5117087 RRF SERVICE DIAGNOSTIC KIT

46-296817P1 50 OHM BNC TERMINATOR

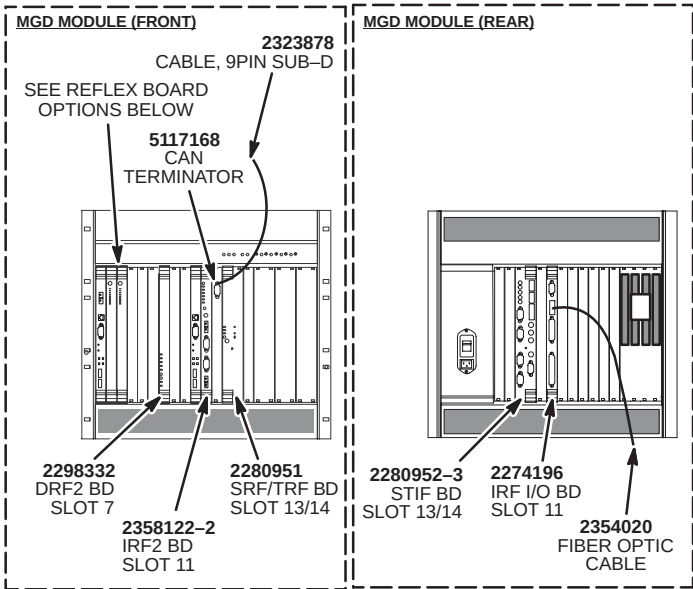

NOTE: The following Array Processor ReFlex board Options are Not Shown.

2325947-12 R400 MINERVA AP (MASTER)

2325947-13 R400 MINERVA AP (SLAVE)

2325947-14 MINERVA PIM

2325947-15 R800 MINERVA AP (MASTER)

2325947-16 R800 MINERVA AP (SLAVE)

NOTE: The following MGD Module parts are Not Shown.

2294300-3 AGP BOARD

2294300-4 APS BOARD WITH IT DAUGHTER BD

2294300-5 BRIDGE BOARD

2294300-6 MGD CHASSIS W/BACKPLANE & AC ENTRY (No Boards)



2294300-7 MGD FAN

2294300-8 or

5181495-3 POWER SUPPLIES

2294300-9 INDICATOR BOARD

2294300-19 SCP2 BOARD

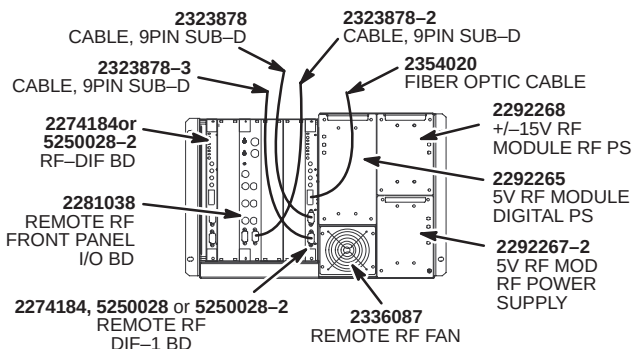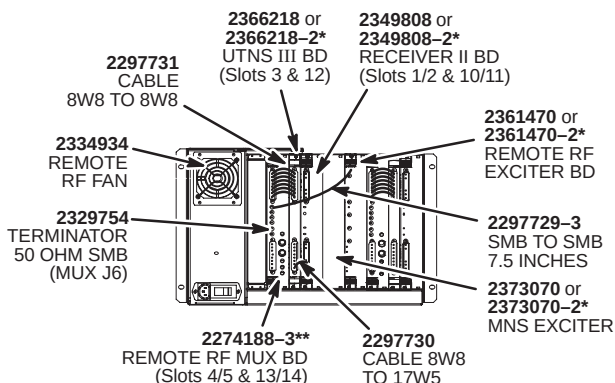2305884 REMOTE RF MOD ASM (CHASSIS W/MIDPLANE) (FRONT)**2305884 REMOTE RF MOD ASM (CHASSIS W/MIDPLANE) (REAR)**

NOTE: Additional REMOTE RF MOD ASM Parts/Cables are located on page 29.

* A part that does not have a "-2" can be replaced with either the same part number or the same number ending with a "-2". If replacing the "-2" version, then only the part number with a "-2" can be used.



** 2274188 and 2274188-2 are replaced by 2274188-3. Also, all 2274188 series MUX boards can be replaced by backwards compatible MUX2 board, 5250106.

NOTE: The following RRF MODULE (page 28) Parts/Cables are Not Shown.

46-306565P1	CABINET FILTER
2294728	COMPACT PCI FRONT PANEL BLANK (8HP) (2 SLOT)
2319152	COMPACT PCI FRONT PANEL BLANK (4HP) (1 SLOT)
2308729	BNC TO SMB CABLE (MUX J8 TO CAB I/F J1)
2308728	8W8 CABLE ASSEMBLY (MUX J5)
2322578	J15/23/25/26 DATA INPUT CABLE (CAB I/F)
2322580	CABLE J18, UNBLANK (CAB I/F)
2322581	CABLE J24, MAGNET MONITOR (CAB I/F)
2322582	CABLE J30/31, SERIAL EXP (CAB I/F)
2322583	CABLE J62 (15 PIN SUB-D) (CAB I/F)
2322586	CABLE J61 (37 PIN SUB-D) (CAB I/F)
2315093	ETHERNET CABLE
2315093-2	CABLE J63 ETHERNET (CAB I/F)
2308729-2	BNC TO SMB CABLE (MUXJ13 TO CAB I/F J2)
2308729-3	BNC TO SMB CABLE (MUXJ14 TO CAB I/F J3)
2394870-3	RRF slot14R-J5 I/O J200
5110796	RRF slot14R-J13 I/O J202
5110795	RRF slot 14R-J13 I/O J202
5115417-3	RRF slot14R-J8 RF Amp J14
5115417	RRF slot12R-J31 Splitter-1
5113193	RRF slot12R-J7 RRF slot11R-J6/7
2394870-2	RRF slot5R-J5 I/O J201
5113193-2	RRF slot3R-J7 RRF slot2R-J6/7
5115417-2	RRF slot3F-J31 Splitter-2
5115630	J7 MUX (2) TO J6 MUX (1) CABLE
5115630-2	RRF slot14R-J11 RRF slot9R-J6


NOTE: The following ASC MODULE (page 30) Parts/Cables are Not Shown.


5115417-4	ASC slot14-J1 RF Amp J8
5115417-5	ASC slot14-J4 RF Amp J7
5115417-6	ASC slot11-J1 RF Amp J6
5120955	ASC POWER SUPPLY
2391148	1.5T BB2 AMP I/F ASM

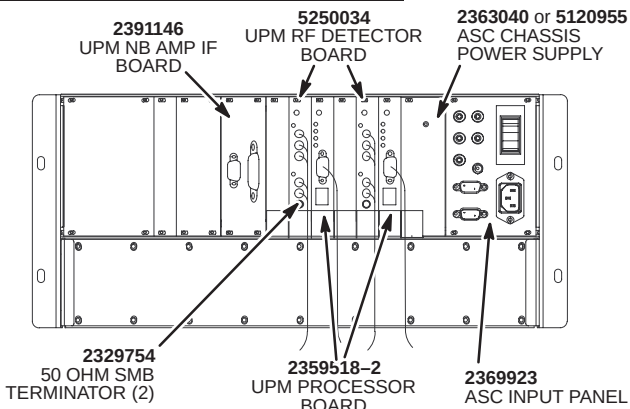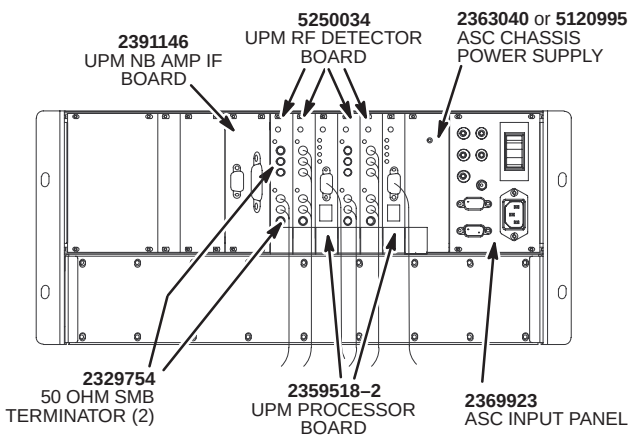2363041 or 2363041-2* ASC CHASSIS ASSEMBLY**2363041 or 2363041-2* ASC CHASSIS ASSEMBLY WITH MNS****NOTES:**


* Additional ASC MODULE Parts/Cables are located on page 29.

[illegible]

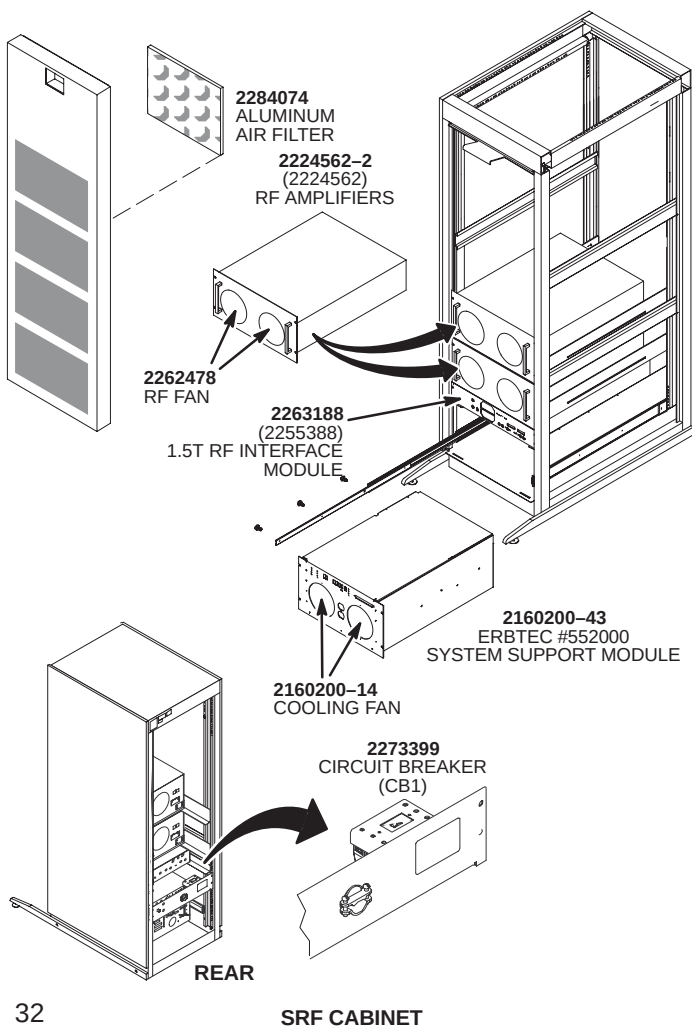EXISTING SRF CABINET UPGRADED TO EXCITE HD

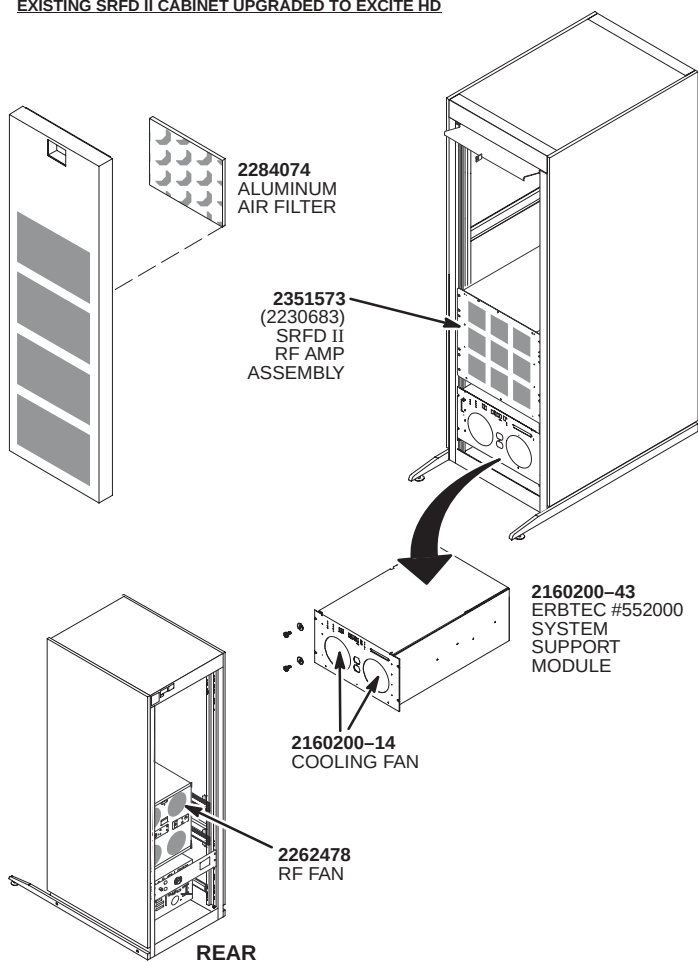
EXISTING SRFD II CABINET UPGRADED TO EXCITE HD


SRFD II CABINET

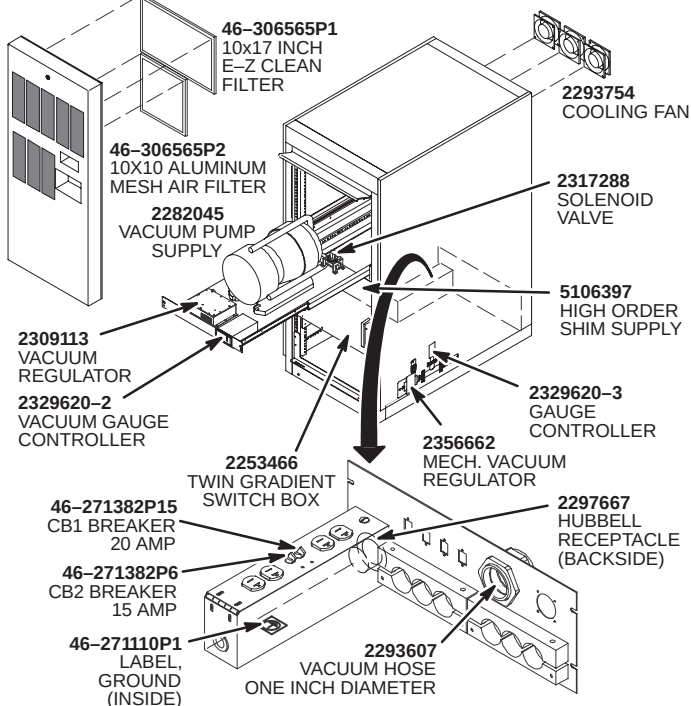2239500 TAC

NOTE: The following parts are Not Shown.

2282045-2	VACUUM PUMP TIP SEAL KIT
2293608	QF25 CENTERING RING & O-RING
2293608-2	QF16 CENTERING RING & O-RING
2293609	QF25 VACUUM ELBOW
2293616	QF25 ADAPTER TO ONE INCH HOSE
2293617	VACUUM REDUCING 'T'
2293618	QF16 CLAMP
2293618-2	QF25 CHAIN CLAMP
2294010-3	HOSE CLAMP, SST ONE INCH
2297760	VACUUM GREASE, MSDS# 8366442
2329620-4	VACUUM GAUGE CONTROLLER CABLE
2284301	VACUUM GAUGE TO TAC J2 CABLE

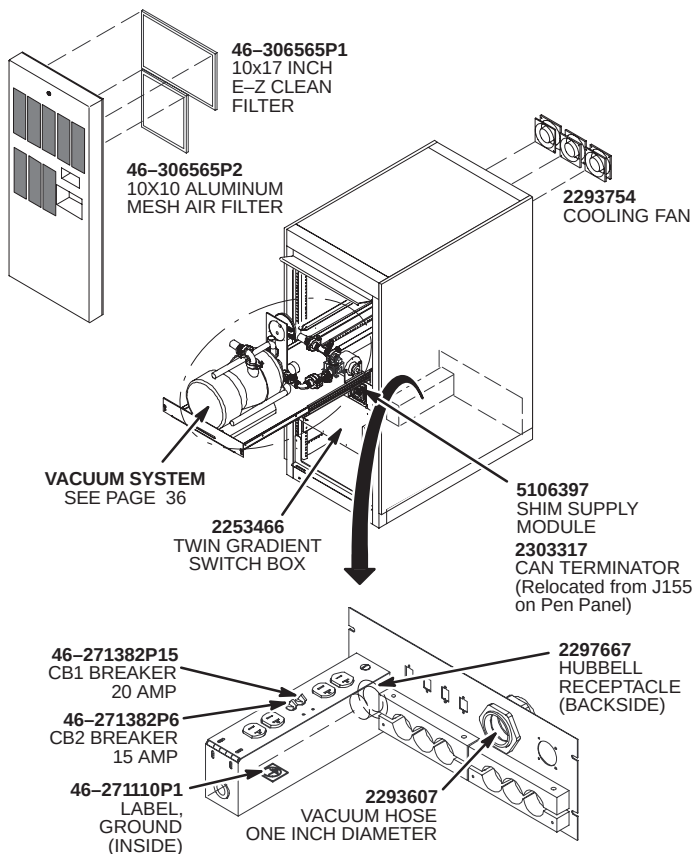
5118162 TAC

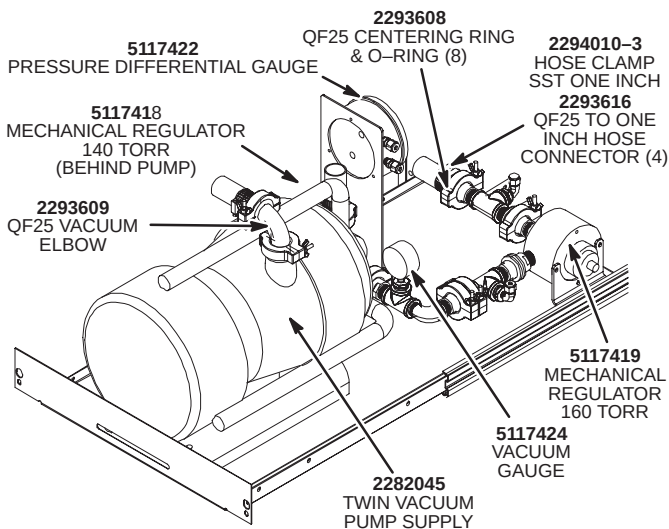

NOTE: The following parts are Not Shown.

2282045-2	VACUUM PUMP TIP SEAL KIT
5118818	CLAMP NW 20-25
2297760	VACUUM GREASE, MSDS# 8366442

5127147
SERVICE KIT FOR TAC VACUUM
SYSTEM UPGRADE
(Used to replace old style vacuum
system, with solenoid)

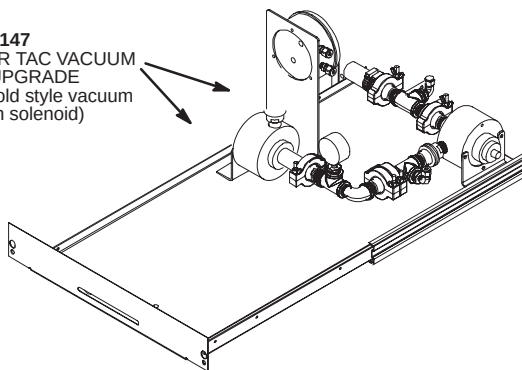

[illegible]

38

NOTES

The table below provides information for generic replacement cables. If any of the cables are ordered, make sure to order write-on labels to be used for identifying cable part number and connection designations at the end of each cable.

GENERIC REPLACEMENT CABLES

Description	Length	
	65 Feet	100 Feet
9 Pin Sub-D Male/Female	46-244714G30	46-244714G37
9 Pin Sub-D Female/Female	46-244713G30	46-244713G37
15 Pin Sub-D Male/Female	46-244712G30	46-244712G37
15 Pin Sub-D Female/Female	46-244711G30	46-244711G37
25 Pin Sub-D Male/Female	46-244710G30	46-244710G37
25 Pin Sub-D Female/Female	46-244709G30	46-244709G37
37 Pin Sub-D Male/Female	46-258877G30	46-258877G37
	50 Feet	65 Feet
BNC RG 223	46-251710G1	46-251710G18
Cable Write-on Labels (60 labels per booklet)	46-208781P3	

DESIGNATED CABLE REPLACEMENTS

CABLE PART #	CABLE RUN
46-244710G30	Run 1227 Run 1236
46-244712G30	Run 1188
46-244713G30	Run 1181
46-244714G30	Run 1180 Run 1181 Run 1233
46-258877G30	Run 1189 Run 1190

CABLE KITS

FRU Kit Description	Part #
Receive Cables for Equipment Room, 34.5 ft (10.56m)	5119953
Receive Cables for Equipment Room, 65.6 ft (20m)	5119953-2
Receive Cables for Magnet Room, 28.2 ft (8.6m)	5119954
Receive Cables for Magnet Room, 49.2 ft (15m)	5119954-2
Transportable/Relocatable MC Rcv Cable Kit (Eqmt. Rm)	5119953-3
Body Receive Cables (both rooms), 34.5 ft (10.52m)	5119956
Body Receive Cables (both rooms), 65.6 ft (20m)	5119956-2

FIBER OPTIC CABLE KIT – FIXED SITE


Description	Part #
Includes clamps, cable ties, mounts, and Runs 1045, 1083, & 1179.	5115713

TIMES MICROWAVE CABLE PARTS

Description	Part #
Straight Male N Connector for LMR 600	2352195
7/8 EIA Right Angle Plug & Center Pin	2319986

ETHERNET CABLE KITS

Description	Part #
Ethernet cables	5160775

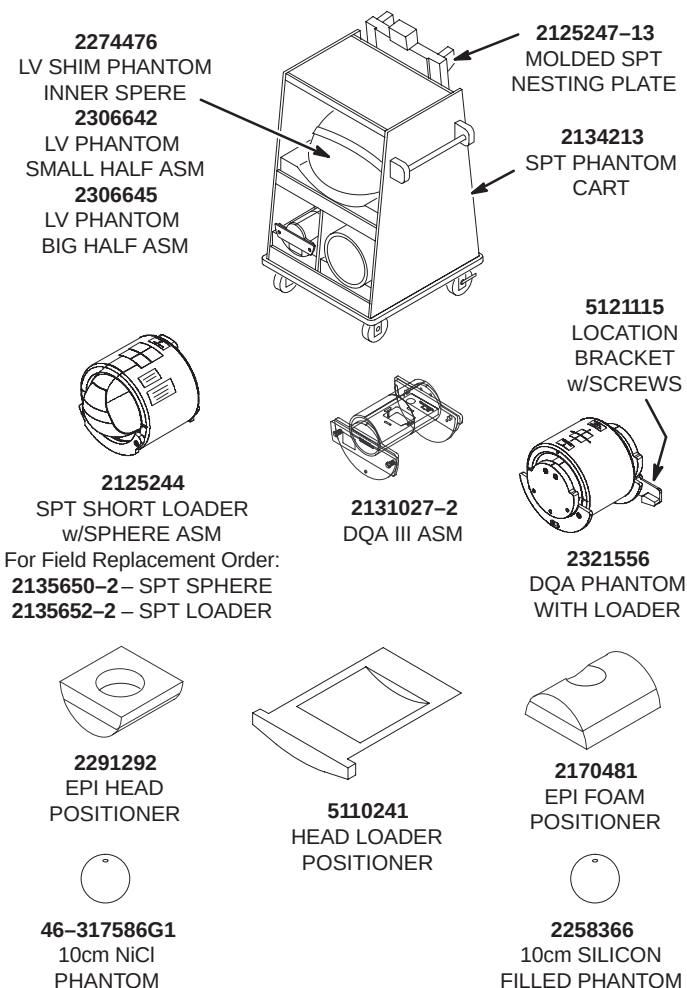

42

NOTES

GRAFIDY KIT:

2270695
CuSO4 SAMPLE VIAL

HEAD COILS:**1.5T SPLIT HEAD COIL ASM**

2341973	1.5T SPLIT HEAD COIL ASM
5191687	HEAD COIL LATCHING FRU KIT
2145391	MIRROR ASM

1.5T GENERAL PURPOSE COIL CABLE W/COIL ID	M1085NJ
2383623	SINGLE CHANNEL RECEIVE ONLY CABLE W/COIL ID

1.5T LINEAR EXTREMITY COIL CABLE W/COIL ID	M1085PC
2383622	LINEAR EXTREMITY COIL CABLE W/COIL ID

1.5T LEGACY HEAD COIL CABLE W/COIL ID	M1085PJ
2384268	LEGACY HEAD COIL CABLE W/COIL ID

EXTREMITY COILS:

1.5T MNS PHOS FLEX COIL	M1040PH
2219090	MNS P31 COIL ASM.
2241998	QUICK DISCONNECT BOX W/CABLE TRAP
46-317287G7	COAXIAL CABLE, 32"

1.5T EXTREMITY SURFACE COIL	M1085AC
46-265442G1	COIL, CYLINDRICAL
46-317278P1	CABLE ASM.

1.5T BODY FLEX COIL	M1085BE
46-328215P1	DUAL-COIL CONNECTOR AND SLEEVE
2151291	ANTERIOR COIL W/CABLE
2151291-2	POSTERIOR COIL W/CABLE

1.5T GENERAL PURPOSE FLEX COIL	M1085GR
2128554	COIL ASM.
2128557	CABLE KIT

1.5T BILATERAL BREAST COIL	M1085BR
46-320785P3	MEDRAD BREAST COIL W/O CABLE
46-320765P2	CABLE TO HEAD PORT
46-320756P2	PA CABLE

1.5T OPEN BREAST COIL	M1087BT
2246360	BREAST COIL
2246360-3	CABLE ASM.

COILS CONTINUED ON NEXT PAGE.

REVISION 7

DIRECTION 5133314

1.5T SHOULDER COIL	M1085AN
46-301058P9	MEDRAD SURFACE COIL
1.5T SHOULDER ARRAY COIL	M1087SD
2225252-7	SMALL COIL ASM.
2225252-6	LARGE COIL ASM.
1.5T SHOULDER PHASED ARRAY COIL	M1087SH
2100937-17	MEDRAD ANTERIOR COIL W/CABLE
2100937-18	MEDRAD ANTERIOR CABLE ASM.
2100937-19	MEDRAD POSTERIOR COIL W/CABLE
2100937-20	MEDRAD POSTERIOR CABLE ASM.
1.5T HIGH RES WRIST ARRAY	M1087HW
2293670-2	HIGH RES WRIST ARRAY W/BASEPLATE ASM.
2293670-3	WRIST PHANTOM
2293670-4	PHANTOM HOLDER
1.5T EXTREMITY COIL W/COIL ID	M1085NA
2383613-2	MAI QUAD KNEE COIL W/CABLE
2383613-3	MAI QUAD KNEE CABLE
1.5T QUAD EXTREMITY KNEE/FOOT COIL	M1085E
2225478-7	VELCRO STRAP/LATCH
2383613-2	COIL ASM.
2383613-3	CABLE ASM.
2225478-10	PHANTOM
1.5T PA KNEE COIL	M1087AA
46-320406P1	COIL ASM.
1.5T PA TR KNEE/FOOT IMAGING OPTION	M1087E
2293674-14	COIL W/CABLE
2293674-13	CABLE
2293674-3	PHANTOM
2293674-5	PHANTOM POSITIONER
1.5T CARDIAC PHASED ARRAY IMAGING OPTION	M1087CA
2174499	CARDIAC ARRAY
2189499	COIL ASM. – ANTERIOR
2189499-2	COIL ASM. – POSTERIOR
2215748	CABLE ASM.
1.5T CARDIAC PERIPHERAL VASCULAR COIL	M1087PR
2293676-2	CVNV PV ARRAY COIL
2293676-3	CABLE ASM.
2293676-4	PHANTOM
2192406-7	LOADING PHANTOM
2293676-6	PHANTOM POSITIONER
2285178-4	PHANTOM STRAP
1.5T PERIPHERAL VASCULAR ARRAY COIL	M1087PW
2192406-2	CABLE ASM
2192406-12	COIL ASM.
2192406-5	PIN DIODES
2192406-6	SNR PHANTOM
2192406-7	LOADING PHANTOM
2192406-9	PHANTOM POSITIONER

COILS CONTINUED ON NEXT PAGE.

REVISION 8

DIRECTION 5133314

1.5T NEUROVASCULAR ARRAY COIL	M1087N
2225476-2	COIL ASM.
2225476-3	CABLE ASM.
2225476-6	COIL LATCHES
2225476-7	PHANTOM POSITIONER
1.5T NEUROVASCULAR ARRAY COIL W/ NEW TOP	M1087NB
2293668-2	COIL W/CABLE
2225476-3	CABLE ASM.
2225476-6	COIL LATCHES
2293668-7	PHANTOM POSITIONER
2293668-8	NV ARRAY CABLE – HD COMPATIBLE
1.5T NEUROVASCULAR HEAD/NECK COIL	M1085HN
2118138-5	COIL ASM.
2383614	QUAD NV CABLE
1.5T ANTERIOR NECK COIL	M1085AK
46-287396P10	MEDIUM SURFACE COIL
1.5T POSTERIOR CERVICAL SPINE COIL	M1085AM
46-287807P12	MEDRAD NECK COIL
1.5T QUADRATURE T/L SPINE COIL	M1085AP
46-307387P10	MEDRAD COIL W/POSITIONER
1.5T QUADRATURE CERVICAL SPINE COIL	M1085AR
46-320422P1	COIL W/CABLE
1.5T CTL COIL	M1087CT
2104799	CTL MULTICOIL ASM.
1.5T CTL SPINE ARRAY COIL	M1087SP
2225545-6	COIL W/CABLE & BRIDGE
2225545-2	CABLE ASM.
2225545-8	ANTERIOR NECK COIL ASM.
2225545-9	ADAPTER BLOCK
1.5T TORSO ARRAY COIL	M1087TB
2215261	COIL
2215748-2	CABLE ASM.
1.5T PHASED ARRAY PELVIC COIL	M1087PA
2144100-2	PELVIC COIL
46-307821G1	CABLE ASM.
2144806	PHASE ADJUST
1.5T TORSO PELVIS COIL	M1087TP
2298043-2	COIL ASM.
2298043-3	CABLE ASM.
2298043-4	PHANTOM
2298043-5	PHANTOM POSITIONER
1.5T HI RES BRAIN COIL	M3087JA
2317112-2	8CHANNEL HR BRAIN CL/BASEPLATE
2317112-3	8CHANNEL HR BRAIN PHANTOM POSITIONER
2317112-8	8CHANNEL HR HEAD CABLE ASSEMBLY
2317112-12	8CHANNEL HR BRAIN LATCH KIT
5146634-6	MIRROR ASM

COILS CONTINUED ON NEXT PAGE.

COILS

REVISION 8

DIRECTION 5133314

1.5T 8 CHANNEL NV ARRAY COIL		M3087JB
2317115-2	8CHANNEL NV ARRAY	
2317115-3	8CHANNEL PHANTOM KIT	
2317115-10	8CHANNEL RIGHT/LEFT LATCH KIT	
2317115-11	8CHANNEL SUPERIOR LATCH KIT	
1.5T 8 CHANNEL CTL COIL		M3087JC
2377426-2	8CHANNEL COIL W/CABLE	
2377426-3	8CHANNEL PHANTOM	
2377426-7	8CHANNEL CABLE	
2377426-10	8 CHANNEL PREAMP	
2418267	ANTERIOR BRIDGE	
2415553	CTL COIL TABLE LOCATOR FRU KIT	
1.5T 8 CHANNEL CARDIAC COIL		M3087JD
2304255-2	8 CHANNEL CARDIAC ARRAY COIL W/CABLE	
2304255-3	8 CHANNEL CARDIAC ARRAY PHANTOM KIT	
2304255-7	8 CHANNEL CARDIAC ARRAY CABLE	
2377426-10	8 CHANNEL PREAMP	
1.5T 8 CHANNEL BODY ARRAY COIL		M3087JE
2377425-2	8 CHANNEL BODY ARRAY COIL W/CABLE	
2377425-3	8 CHANNEL BODY ARRAY PHANTOM KIT	
2377425-7	8 CHANNEL BODY ARRAY CABLE	
2377426-10	8 CHANNEL PREAMP	
1.5T 8 CHANNEL T/R KNEE COIL		M3087JF
5114356-2	8 CHANNEL T/R KNEE ARRAY COIL W/CABLE	
5114356-3	8 CHANNEL T/R KNEE ARRAY CABLE	
5114356-4	8 CHANNEL T/R KNEE ARRAY BASEPLATE	
5114356-11	8 CHANNEL T/R KNEE ARRAY PHANTOM	
5114356-12	8 CHANNEL T/R KNEE ARRAY PHANTOM	
5114356-13	8 CHANNEL T/R KNEE ARRAY PHANTOM POSITIONER	
1.5T BILATERAL BREAST COIL		M3087JG
2412286	BREAST COIL W/O CABLE	
2412283	BREAST COIL CABLE	
2412281	BREAST COIL HEADREST ASM	
2412282	BREAST COIL QA PHANTOM	
2412285	BREAST COIL BIOPSY ARM	
2414541	BREAST COIL BIOPSY MOUNT	
2408500	BREAST COIL PATIENT PAD	
2408503	BREAST COIL CHAMBER COVER	
1.5T HD NV ARRAY COIL		M3087JJ
51107092-2	HD NV ARRAY COIL W/CABLE	
51107092-3	HD NV ARRAY CABLE	
51107092-6	HD NV PHANTOM & PHANTOM POSITIONER	
51107092-7	HD NV ARRAY MIRROR REPLACEMENT	

[illegible]

MNS OPTIONS:

Description	Part #
1.5T 31P TR Switch	2354050-2
1.5T 13C TR Switch	2354050-3
1.5T 13C MNS Preamp	5114799
MNS QA Probe 1.5T 13C	5114403
MNS QA Probe 1.5T 31P	5114978-2
1.5T 13C Sample Vial	5114978-3
1.5T 31P Sample Vial	5114978-4

OPTIONS:

CONFIRMA MAMMO CAD M3033NF
5117875 COMPLETE PC WITH SOFTWARE

EXCITE HD TWIN RESISTIVE SHIM POWER SUPPLY M3000PP
5118323 1.5T HIGHER ORDER SHIM POWER SUPPLY KIT

CLARIVIEW IMAGE FILTERING M1090PN
2203037 AW1.2 CLARIVIEW FLOPPY ASM.
2203038 AW2.0 CLARIVIEW FLOPPY ASM.
2203039 SECURITY KEY FLOPPY ASM.

BAR CODE READER KIT
2180690-2 READER_ONLY
2180690-3 Y CABLE_ONLY

ReportCARD HARDWARE
5117866-26 ReportCARD COMPUTER
2390367-4 ReportCARD SOFTWARE COLLECTOR
2390369 ReportCARD OPTON KEY

EXCITE SOFTWARE OPTION KEYS:

2156158-60	Release 12.x Software Collector EXCITE HD (M4)
2383494	Fiesta C
2383495	3D Delayed Enhancement
2383496	Navigator 3D
2383497	3D FatSat Fiesta
2383498	Tricks
2390037	Asset
2390040	ConnectPro
2390075	3D FRFSE
2390078	EchoPlus
2390079	Diffusion Tensor
2390084	Probe 2000
2390085	Probe 200 Upgrade
2390086	Sage 7
2390091	2D Delayed Enhancement
2390095	iDrive ProPlus
2390096	Cardiac Tagging
2390097	SmartPrep 2000
2390098	Flouro Trigger MRA
2390103	Performed Proc Step
5110228	Single Voxel Probe-P
5110229	Single Voxel Probe-S and Probe-P
5110230	Probe 2D – CSI
5110231	Probe 3D – CSI
5110232	MNS, SAGE, Probe-P
5110319	Research Proton Spectro
5120544	Release 12.0 M4 Apps DVD
5121548	ReportCard on AW
5121550	Tractography on AW
5122363	LAVA
5122364	PFO Evaluation
5122365	PROPELLOR Brain Suite
5122366	Tractography
5122367	2D FatSat FIESTA
5122368	3D Flair
5122370	RealTime FIESTA
5122371	VIBRANT-II
5122372	8 Channel EXCITE 3
5122373	16 Channel EXCITE 3
5122374	1.5T EXCITE ScanTools 12.0

SOFTWARE OPTIONS

BrainWave OPTIONS:

2363727	15" LCD Display with 5 foot VGA Cable for Computer connection
2394034-5	(Cedrus) Response Pad Asm. Left & right hand units with two buttons for each hand. Includes fiber optic cable (20 ft) and optical converter box
2394034-6	(Cedrus) "Lumina" response controller box
2394034-7	(Cedrus) AC Power module for Lumina Controller Box
2394034-9	FRU usage only Kit containing all cables [Cedrus-4, -5, -6, -7, 2xxxxxx-6, 46-265505P1 (qty 2)]
2394034-10	Penetration Sub-Panel Asm
2394034-16	PC compatible Keyboard (narrow)
2394034-17	PC compatible Mouse
5110767	BrainWave Application Software CD-ROM
5110736	EXCITE II Scanner License Key MOD for the BrainWave H/W Option
5110946	PC and Platform Software Asm
5110968	Patch CD for BrainWave Stimulus PC

[illegible]

REVISION HISTORY

REV/DATE REASON FOR CHANGE

- 1 Nov 16, 04 Initial release
- 2 May 4, 05 Corrected RF coil, temp sensor cable, and split bridge part numbers. Added SRI-3 upgrade kit. Added CAN terminator to PHPS. Corrected Host monitor and SCSI tower part numbers. Updated HFD and HFD-S/PDU cabinet modules. Added RF Amp FRU. Added SRF and SRFD2 cabinets. Added CAN terminator and vacuum upgrade kit to TAC cabinet. Updated MDP and cable replacement information.
- 3 Nov 16, 05 Corrected front split bridge asm part number & isolation box description. Added new part numbers for small & large cable tracks for long bore magnets. For OW, added DVD upgrade kit, video cable, and notes for HP800 & HP8200 substitution rules. Updated RF Mux board part number. Added HOS Supply to TAC. Added fiber optic cable kit. phantoms & new loader parts, split head coil, and NV HD array cable. Corrected and updated other parts as necessary.
- 4 Aug 1, 06 Corrected part #'s for gradient clamp block, 8W8 cable in Rear Ped & LPCA cover. Added FRU's to Patient Table. Added ethernet hub to GOC and ethernet cable kit. Added HD NV Array Coil. Added Cardiac Tagging s/w option key (MRIhc12672).
- 5 Jan 26, 07 Added RF Shield Change Out Kit, Speaker Grill nylon screw, LPCA Wheel Housing FRU kit (MRIhc21559), & Dock pedal FRU parts (MRIhc17629). Added new LCD Host Monitor, ASC Module UPM board and RRF Module MUX2 board. Removed MCR Tool.
- 6 Mar 13, 07 Added front cover laser, 5 gallon coolant pail, and LPCA LED cable. Added note for referring to 2138247 for parts not in this manual. Added center pin to description of Times Microwave right angle connector (CQA 13099582). Removed ECMT kit (CQA 13097585). Added Head Coil Latch Kit.
- 7 Sep 17, 07 Moved Revision History to last page. Added new Netcom Hybrid Splitter and FRU's, revised existing Teledyne splitter FRU's. Updated HFA-S Asm part # in HFD-S cabinet (ECO 20428110). Revised MGD & ASC Module part #'s. Added BRM gray insulation foam (PQR 13047306), LPCA bracket, pillar, and screws FRU kit (PQR 13121160) & updated LPCA cover part # (PQR 13083953), added Patient Table bumper and panel FRU kit (PQR 13127363), Grafidy CuSO4 vial (PQR 13098795), and 8-Channel CTL anterior bridge (PQR 13115149).
- 8 May 28, 08 Added Ground Strap kit part #. Rear Pedestal: Updated Longitudinal Encoder part number (PQR 13164271) and added Encoder board part # 46-226552G2. Chiller: Added E & W mobile Heat Exchanger part #'s (PQR13133952). Added Lytron pump coupler part # (PQR 13136000). HFD: Updated FHA Control Bd part #. MGD: Updated SCP2 Bd part #. Added 5181495-3 MGD P.S. Coils: Added Mirror assembly to Brain Coil (ECR 203882). Added M3087JG Bilateral Breast Coil. Option Keys: Updated EXCITE HD (M4) software part #.
- 9 Jun 15, 09 Added Flexible Coupling (PQR 13206931) to Rear Pedestal. Updated Monitor offerings on OW (PCN 290188). Also added HP8200 DIMM (ECR 2078720) and updated to 73 Gb HDD (ECR 2065991). In RRF Chassis, added 5250028-2 RF-DIF (ECO 2080729).
- 10 Nov 3, 09 Added LPCA Multicoil Asm, 5400062, and Preamp Bd, 5400044 (ECO 2083575). Added new GOC AA 2395043-2 (ECO 2084000).