

Operator's Manual

1.5T HD Open Breast Array for the GE Value MR System



Manufacturer
Invivo
Invivo Corporation
3545 SW 47th Avenue
Gainesville, FL 32608 U.S.A.
500166 Rev 07

CE 0123

Distributed by:
GE Healthcare



GE Part Number 5145567-3

First Edition
July 2005
©2005 Invivo Corporation

All rights reserved. No part of this publication may be reproduced, transmitted, transcribed, stored in a retrieval system, or translated into any language in any form by any means without the written permission of Invivo Corporation.

Licenses and Trademarks

The Invivo Logo is a registered trademark of Invivo Corporation. "Value" is a registered trademark of the General Electric Company.

Proper performance of this coil is guaranteed only while the coil is being used on the MR system (hardware/software level) specified at the time of purchase. Upgrades or other modifications to the system software and/or hardware may affect compatibility. Prior to upgrading your MR system, please contact your GE Healthcare representative to discuss coil compatibility issues. Failure to do so may void your warranty.

Medical Device Directive

These products conform with the requirements of council directive 93/42/EEC concerning medical devices, when they bear the following CE Mark of Conformity:

European Representative

GBM Authorised Representative Ltd.
The White House, 2 Meadow
Godalming, Surrey GU7 3HN
United Kingdom



 Attention, Consult Accompanying Documents	NOTICE: This equipment shall be transported and stored under the following conditions: <u>Shipping</u> 1. Ambient temperature range of +5°C to +65°C. 2. Relative humidity range of 10% to 100%, non-condensing. <u>Storage</u> 1. Ambient temperature range of +10°C to +40°C. 2. Relative humidity range of 10% to 90%, non-condensing. Atmospheric pressure range of 500 hPa to 1060 hPa.
 Type BF Equipment	
 Class II Equipment	

WARNING: This product contains chemicals, including lead, known to the state of California to cause birth defects or other reproductive harm. **Wash hands after handling.**

Caution:
Federal law restricts this device to sale, distribution and use by or on the order of a physician.

INTRODUCTION

This manual describes the safety precautions, features, use and care of the Invivo 1.5T HD Open Breast Array, for use with the GE Value MR System. Please review this manual thoroughly before using the device.

If you have any questions or comments regarding this manual, or need any assistance with the use of the product, please contact your GE Healthcare representative at:

1-800-GE CARES

COMPATIBILITY

The Invivo 1.5T HD Open Breast Array, Invivo model number 800358, is compatible with the GE Value MR System. The minimum required system software level is 14.0.

EXPLANATION OF SYMBOLS



Medical Device Directive - These products conform with the requirements of council directive 93/42/EEC concerning medical devices, when they bear the CE Mark of Conformity.



Attention, Consult Accompanying Documents.



Type BF Equipment




Class II Equipment



Do not cross or loop cables. Arcing and patient burns could result.



Waste in Electrical and Electronic Equipment

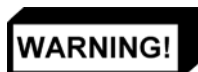

A hazard exists which could cause serious injury.




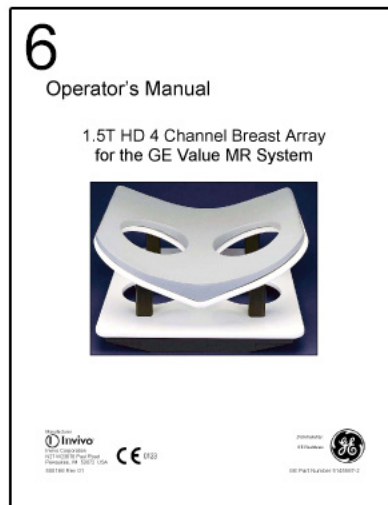


TABLE OF CONTENTS

INTRODUCTION	3
COMPATIBILITY	3
EXPLANATION OF SYMBOLS	3
CHAPTER 1 - 1.5T HD OPEN BREAST ARRAY CONTENTS	5
CHAPTER 2 - SAFETY	6
2-1 Training	6
2-2 Quality Assurance	6
2-3 Indications	6
2-4 Contraindications	6
2-5 Precautions	6
2-6 Cautions	7
2-7 Emergency Procedures	8
2-8 Technical Considerations	8
CHAPTER 3 – QUALITY ASSURANCE	9
3-1 Tools Required	9
3-2 Setup	9
3-3 Procedure	11
CHAPTER 4 - USING THE 1.5T HD OPEN BREAST ARRAY	12
4-1 Positioning the 1.5T HD Open Breast Array	12
4-2 Ramp and Sternum Pad	13
4-3 Positioning the Patient	13
4-4 Connecting the Cable	14
4-5 Establishing the Patient Landmark	14
4-6 Field of View and Coverage	14
4-7 Patient Hearing Protection	14
CHAPTER 5 – SCANNING	15
5-1 Selecting the Coil	15
5-2 Auto Pre-Scan	15
5-3 Localizing with the Body Coil	15
5-4 Using AUTOSHIM	15
5-5 Fat Saturation	16
5-6 Scan Protocols	16
CHAPTER 6 – INSTALLATION AND MAINTENANCE	17
6-1 Installation and Configuration	17
6-2 Cleaning	17
6-3 Storage	17
6-4 Replaceable Accessories	17
6-5 Contact Information	18
CHAPTER 7 – ADDENDUM	19

CHAPTER 1 - 1.5T HD OPEN BREAST ARRAY CONTENTS

The 1.5T HD Open Breast Array consists of the following parts. Please inspect upon receipt to make sure all parts have arrived and are in good order. Use this guide to refer to part names throughout this manual.

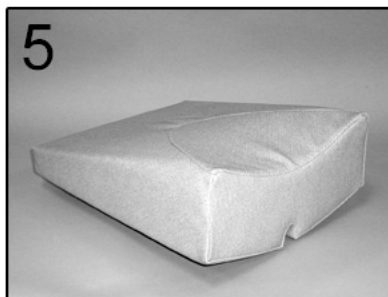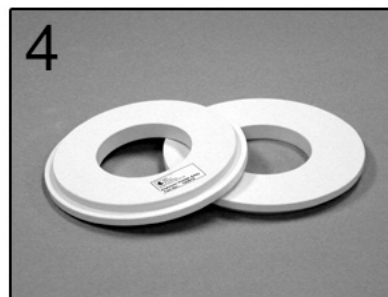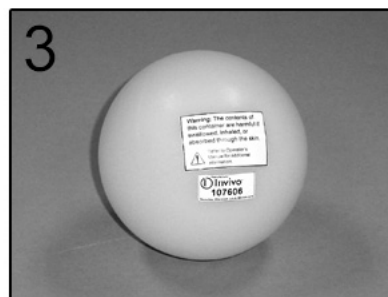

- | | | |
|-----------------|----------------------------|----------------------|
| 1. Breast Array | 3. 4.38" Spherical Phantom | 5. Ramp Pad |
| 2. Sternum Pad | 4. Phantom Positioner (2) | 6. Operator's Manual |

CHAPTER 2 - SAFETY

2-1 Training

This manual contains detailed information regarding the set-up, positioning and use of the Invivo Corporation coil. The instructions must be read carefully and thoroughly before attempting to scan patients with this coil.

2-2 Quality Assurance

The procedure described in the Quality Assurance chapter of this manual should be performed upon receipt of the coil to establish a baseline of coil performance. The procedure should be repeated at regular intervals.

2-3 Indications

The coil is indicated for use, on the order of a physician, in conjunction with a MR scanner as an accessory to produce images of the breast.

2-4 Contraindications

The operator should be aware of the following contraindications for use related to the strong magnetic field of the MR system:

 Scanning is contraindicated for patients who have electrically, magnetically or mechanically activated implants (cardiac pacemakers, for example). The magnetic and electromagnetic fields produced by the MR system and coil may interfere with the operations of these devices.

 Scanning patients with intracranial aneurysm clips is contraindicated.

2-5 Precautions

Precautions should be taken when scanning patients with the following conditions:

 Greater than normal potential for cardiac arrest

 An increased likelihood for developing seizures or claustrophobia

 Unconscious, heavily sedated, or confused patients

 Inability to maintain reliable communications

2-6 Cautions

The following general warning statements apply to scanning with a magnetic resonance system. For further details, review the warnings in your MR system Operator's Manual.

  Do not cross or loop cables. Arcing and patient burns could result. Route cables out of the magnet so that they do not touch the patient.

 Ensure that the patient does not form a loop with any body parts. For breast coils in particular, patients are imaged lying prone with both arms extended above the head (superman). **Do not allow patients to touch any part of the right hand or arm to any part of the left hand or arm. The loop formed by doing so could cause an RF burn at the point of contact. Use pads to ensure the patient's hands do not touch any metal features of the patient table or cradle.** Be sure to educate the patient accordingly and check the patient's position immediately before the scan.



 Assure that the patient is not touching the bore. If necessary, place pads between the patient and the surface of the bore.

 If the patient complains of warming, tingling, stinging, or similar sensations, promptly stop the scan procedure, examine the patient, and contact the responsible physician before continuing the procedure. Pay special attention to very young, sedated, or other compromised patients who may not be able to communicate effectively.

 Patients with implanted ferromagnetic metal should not be scanned because the magnetic field may interact with implanted surgical clips or other ferromagnetic materials.

 Persons with cardiac pacemakers or other implanted electronic devices should not enter the magnetic field zone delineated by the MR system manufacturer.

 There is a risk to scanning feverish or decompensated cardiac patients.

 Facial makeup should be removed before scanning because it may contain metal flakes which can cause skin and eye irritation. Permanent eyeliner tattoos may cause eye irritation due to ferromagnetic particles.

 Patients who work in environments in which there is a risk of having embedded metallic fragments in or near the eye should be carefully screened before undergoing a MR exam.

 Visually inspect the cable insulator jackets, strain reliefs, and connector boxes before each use. If the insulation is broken, or if the cable is frayed, immediately discontinue use of the device.

2-7 Emergency Procedures

In the unlikely event that a coil creates smoke, sparks or makes an unusually loud noise, or if the patient requires emergency assistance:

- Stop the scan if one is in progress.
- Remove the patient from the scan room if medical treatment is needed.

2-8 Technical Considerations

-  The coil and accessories require special conditions regarding electromagnetic compatibility. The coil must be installed and used in a shielded scan room provided with the MR magnet and system. The user must ensure that the scan room door is closed during system use. Failure to do so may cause reciprocal interference with portable and mobile RF communications equipment, affecting the performance of the MR coil and/or such equipment.
-  The coil should only be used with the accessories specified in the Operator's Manual.
-  The use of accessories other than those specified in the Operator's Manual may result in decreased ESD immunity of the coil or MR system, causing damage to the coil and/or system.
-  The equipment should not be used with other coils or equipment present in the MR scanner except as specified in the Operator's Manual.
-  Tampering with the cable pins and connector may damage the connector and affect coil or system performance. Please verify that connector and pins are not damaged before use.

CHAPTER 3 – QUALITY ASSURANCE

3-1 Tools Required


TOOLS REQUIRED – TABLE 3-2-1

Description	GE Part #	Invivo #	Qty
4.38" Spherical Phantom	46-317586G1	N/A	1*
4.38" Spherical Phantom	46-317586G1	107606	1*
Phantom Positioner	2246360-9	100679	1

* One phantom is provided with the Value system, and one is provided with the Breast Array package.

3-2 Setup


Place the 1.5T HD Open Breast Array on the patient table with the cable toward the magnet. Note the serial number of the coil.

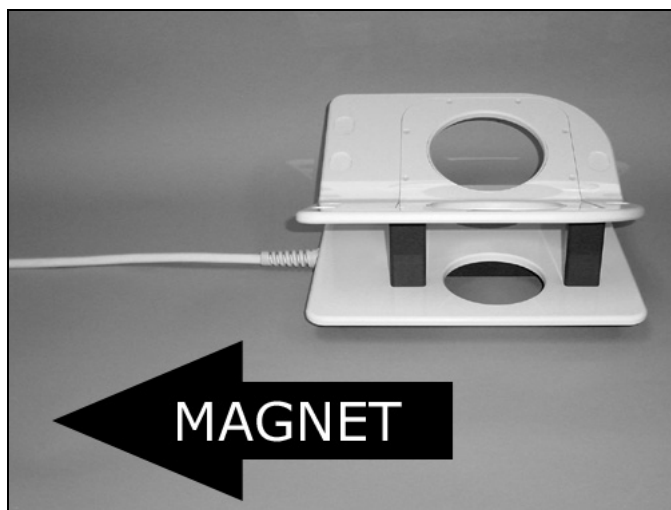

Connect the 1.5T HD Open Breast Array to the system using the Port A connector.




Place one phantom positioner in the right aperture, and one in the left aperture. The phantom positioner must be inserted into the aperture with the flat side facing up.

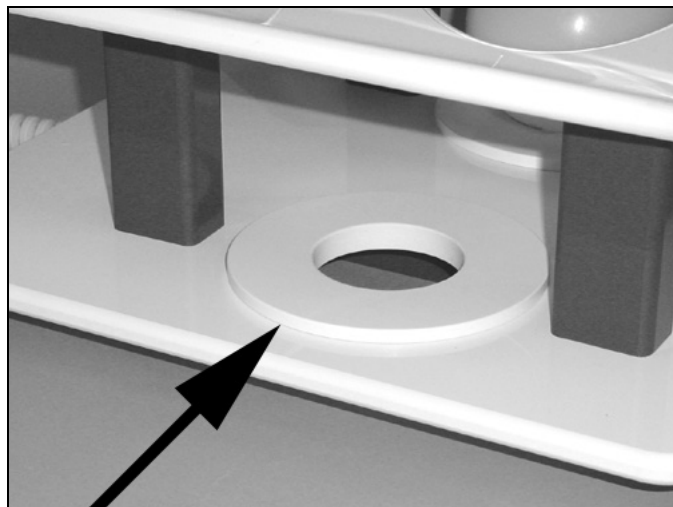

Place one 4.38" spherical phantom on the right phantom holder and one on the left phantom holder.



Establish an axial landmark on the center of the phantoms. Advance to scan.

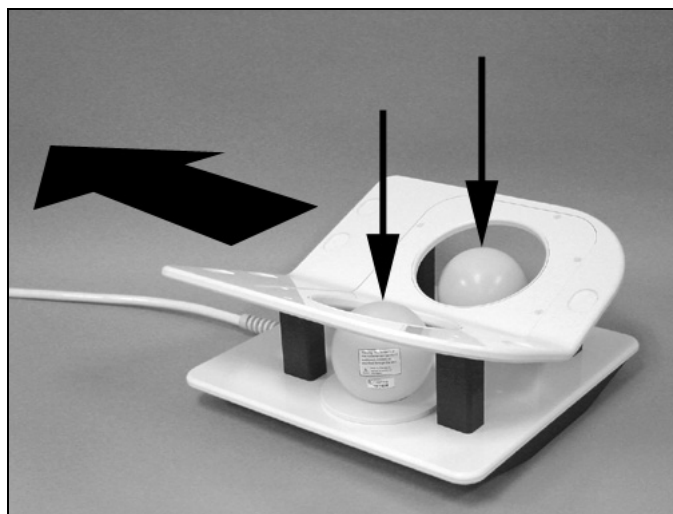

3-3 Procedure


Be sure that you have chosen the correct coil, and activate the desired coil(s) from the menu selection on the surface coil screen. The coil(s) chosen must be R_BREAST, L_BREAST, or 2_BREAST. Activate one of the above coil settings from the operator's console.

Use the following spin echo sequence, TE=20 mS., TR=300 mS., 32 cm. FOV axial, 5 mm. thick, 5 mm. spacing, 256x128, and one average or equivalent sequence. Scan using Auto Prescan. If Auto Prescan does not complete successfully, do not use the coil clinically. Upon completion of Prescan, execute the scan. Window the scan as you normally would (window levels vary from system to system depending upon your hardware revision level). The scan should appear uniform.

If the image has signal voids, distortions, or black streaks, discontinue use of the coil. If the image is excessively noisy, repeat the scan using the GE 5-inch general purpose coil (or alternate coil) supplied with your system. If the 5-inch coil image exhibits similar artifacts, the problem likely lies with the MRI system.

Using the 4.38-inch diameter sphere phantom, scan using the protocol shown.

If the image appears as shown, the coil is malfunctioning. Discontinue use of the coil and contact your GE Healthcare representative.

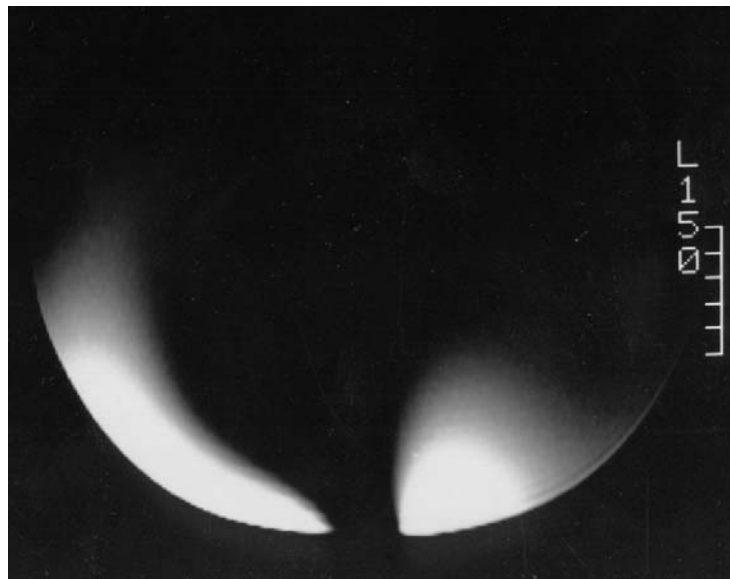

Mode : Multi
PSeq: ME


TR: 300
TE: 12.0 1/1

256x128/1.0 NEX

FOV: 14 cm.
Thk: 10.0 mm .
Imgs: 5/00:42

CHAPTER 4 - USING THE 1.5T HD OPEN BREAST ARRAY

4-1 Positioning the 1.5T HD Open Breast Array

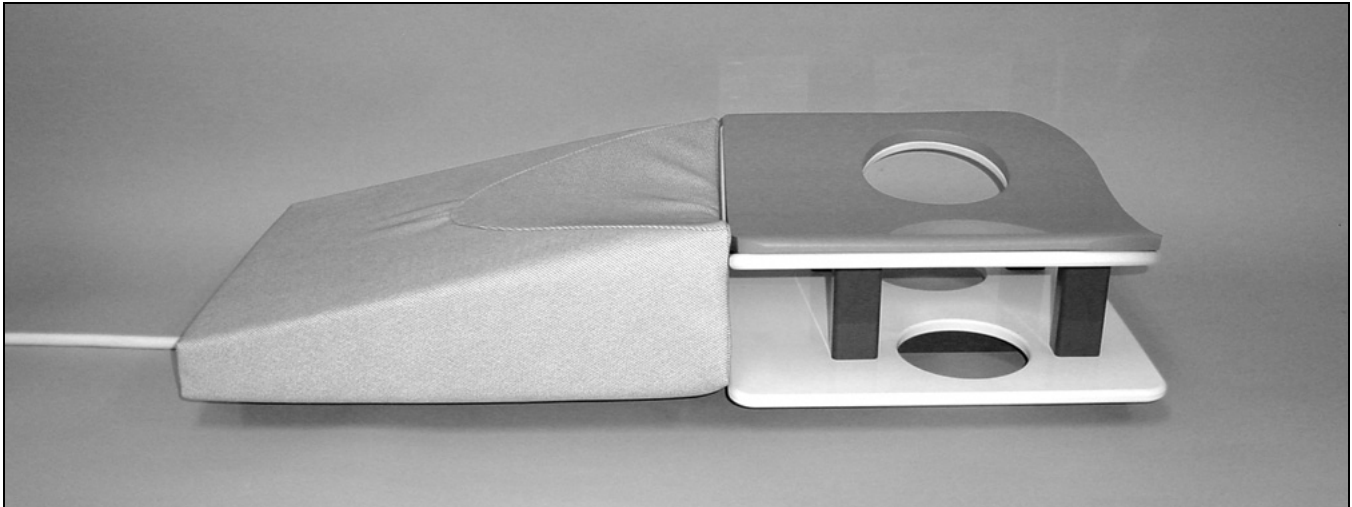

The 1.5T HD Open Breast Array is designed for unilateral imaging of either the left or the right breast, or bilateral imaging of both breasts.

The coil is designed to rest directly on the patient table. It is important that the coil lies directly on the patient table, since this will elevate the patient's breast to the iso center of the coronal scan plane. The Value system, with 14.0 software, allows a scan offset in only one direction. Since the breast is already offset in the lateral direction, it is important that the anterior/posterior (A/P) location of the coil is at zero (in the coronal plane).


If a foreign object (sheets, straps, etc.) is beneath the coil housing, the coil may move during the image acquisition, resulting in motion artifact. The coil must be placed on the patient table with the cable toward the magnet.

- Place the 1.5T HD Open Breast Array on the patient table with the cable toward the magnet.
- The 1.5T HD Open Breast Array is designed for feet first exams only.
- The 1.5T HD Open Breast Array must be placed near the handle end of the patient table to allow the patient to enter the scanner feet first.

4-2 Ramp and Sternum Pad

For patient comfort, a ramp pad and sternum pad have been provided. The large end of the ramp pad (1) is positioned flush with the cable end of the breast platform. The cable is positioned into the groove (2) underneath the ramp pad. The small end of the ramp pad (3) fits flush with GE table pads.

The sternum pad attaches to the top of the breast platform (4) using velcro fasteners. Center the sternum pad openings over the apertures of the breast platform and apply the pad to the platform.

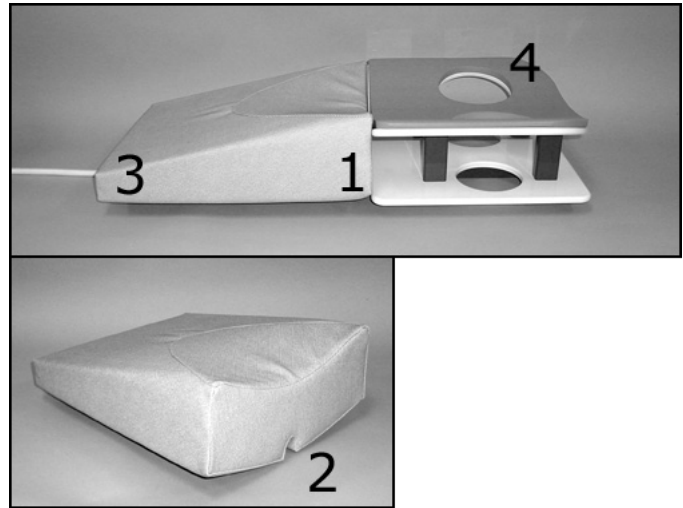

If the ramp or sternum pads become soiled, they may be cleaned by using a solution of 30% isopropyl alcohol and 70% tap water.

4-3 Positioning the Patient

Position the patient prone, feet-first, centering the breasts in the apertures of the coil. Use the open lateral side of the coil to visually center the breasts in the coil.



Ensure that the patient does not form a loop with any body parts. For breast coils in particular, patients are imaged lying prone with both arms extended above the head (superman). **Do not allow patients to touch any part of the right hand or arm to any part of the left hand or arm. The loop formed by doing so could cause an RF burn at the point of contact. Use pads to ensure the patient's hands do not touch any metal features of the patient table or cradle.** Be sure to educate the patient accordingly and check the patient's position immediately before the scan.



4-4 Connecting the Cable

Connect the 1.5T HD Open Breast Array to the system using the standard phased array connector.



4-5 Establishing the Patient Landmark

Use the open lateral side of the coil to visually center the breast in the coil. Position the patient slightly superior or inferior until the breast is centered in the aperture. The patient must be positioned feet first.

4-6 Field of View and Coverage

For unilateral imaging, a FOV of 16 cm. to 20 cm. is suggested, depending upon breast size. The 1.5T HD Open Breast Array aperture is 160 millimeters (16 cm.) in diameter. If both coils are active, a field of view of 32 cm. to 36 cm. is suggested. You may use a smaller FOV if desired. It is recommended that the NO PHASE WRAP option be selected when establishing the scan parameters.

4-7 Patient Hearing Protection

Provide ear plugs for the patient after all instructions have been given.

WARNING!

Hearing protection is required for all people in the scan room during a scan to prevent hearing impairment. Acoustic levels may exceed 99 dB(A).

Hearing protection must have a Noise Reduction Rating (NRR) of 28 dB or better (e.g., 30 dB, 32 dB, etc.).

CHAPTER 5 – SCANNING

5-1 Selecting the Coil

The coil description displayed to the user is HD 4Ch Breast Array by Invivo. The breast imaging platform contains both a left breast coil and a right breast coil. One or both coil(s) can be active at any given time. Make the proper selection from the coil selection screen.

5-2 Auto Pre-Scan

Generally, the Value System auto pre-scan will function normally when used with the 1.5T HD Open Breast Array. If auto pre-scan returns the message, "auto pre-scan failed to complete normally," then manual pre-scan may be required to acquire the image. Consult your GE Value Operator's Manual for manual pre-scan procedures.

If a "TR DRIVER FAULT" occurs, reseal the cable connector and try again. If the fault recurs, then it is likely that the wrong coil selection was chosen at the coil menu. If the problem persists, call your GE Service Representative.

5-3 Localizing with the Body Coil

The 1.5T HD Open Breast Array is designed to allow the Value body coil to be used while the 1.5T HD Open Breast Array is in the scanner. This allows a large FOV (40 cm.) body coil localizer, which is helpful in determining the left or right offset required for imaging the breast.

5-4 Using AUTOSHIM

AUTOSHIM improves image quality by improving the magnetic field homogeneity within the FOV selected. The improvement in image quality is often dramatic when the selected FOV is far off center, and when acquiring FSE or T2 weighted scans.

It is recommended that auto shim be selected only when operating this coil in the unilateral mode. For use in the bilateral mode, turn AUTOSHIM off. To use AUTOSHIM with the 1.5T HD Open Breast Array in the unilateral mode, use the following procedure:

1. Begin a new exam.
2. Choose L_BREAST or R_BREAST as the coil selection after it has been programmed on the system. This is a 1 channel selection.
3. Make sure AUTOSHIM is selected.
4. AUTOPRESCAN.
5. Carry out a series of images for scout purposes, if desired. All subsequent scans will use these AUTOSHIM settings.
6. When switching to the other breast, repeat steps 2-5.

5-5 Fat Saturation

Unilateral imaging is preferred when using FAT SAT, WATER SAT, or SILICONE SAT, since the saturation pulse sequences generally will not work well when the FOV encompasses both breasts due to body coil inhomogeneities. By programming a coil selection key for left breast and right breast, using the appropriate channel for applying a saturation pulse during manual or auto pre-scan is assured.

5-6 Scan Protocols

Invivo Corporation recommends that you select imaging protocols that have been established by your radiologist. Additional protocols may be found within your GE Application Guides.

CHAPTER 6 – INSTALLATION AND MAINTENANCE

6-1 Installation and Configuration

Coil installation and configuration requires the services of your GE Service Representative. Please refer to the Service Manual for coil installation and configuration information.

6-2 Cleaning

The 1.5T HD Open Breast Array, base and pads must be cleaned and stored using the following procedure:

- Wipe with a cloth that has been dampened in a solution of 10% bleach and 90% tap water, or 30% isopropyl alcohol and 70% tap water.
- **Do not pour any cleaning solution directly on the coil!**
- Let the coil housing and pads dry before use.
- **Under no circumstances should the coil be placed into any type of sterilizer.**

6-3 Storage

This equipment must be stored under the following conditions:

- Ambient temperature range of +10°C to +40°C
- Relative humidity range of 10% To 90%, non-condensing
- Atmospheric pressure range of 500 hPa TO 1060 hPa

To store the coil, a storage space of 42.4 cm / 16.7 in (width) x 40.6 cm / 16.0 in (depth) x 18.6 cm / 7.3 in (height) is required.

6-4 Replaceable Accessories

The following table lists replaceable accessories that may be purchased from Invivo Corporation:

Description	Invivo #
Ramp Pad	870254
Sternum Pad	870255

6-5 Contact Information

For additional information, contact:

GE Healthcare Americas
USA 800-558-5102
Canada 800-668-0732

GE Healthcare Asia
China 86-21-62192228
Taiwan 886-2-2505-7900
Singapore 65-291-8528
Australia 61-2-9975-5501
Japan 81-120-48-2630
Korea 82-31-740-6119
India 91-80-845-2923

GE Healthcare Europe
(33) 1-41-19-76-76

GE Healthcare Web Site
www.gehealthcare.com

Invivo Corporation
800-524-1476
262-524-1402

Invivo Web Site
www.invivocorp.com

CHAPTER 7 – ADDENDUM



注释：可替换的识别标签已随模提供，标签规定了模的使用方法，并附有多种语言版本。请依实际情况选择合适的标签直接贴于已有标签之上。

关于模的多语言模标签：106764

BEMÆRK: Skiftende identifikationsmærkater med håndteringsvejledning på flere sprog, leveres med dette fantomaggregat. Vælg den passende mærkat og placer den direkte ovenpå den eksisterende mærkat.

Flersproget mærkatark til fantom: 106764

N.B.: Bij deze fantoomassemblage zijn wisselende identificatielabels geleverd die in een verscheidenheid aan talen instructies geven voor hantering van het fantoom. Selecteer het juiste label en plaats dit direct bovenop het bestaande label.

Meertalig labelsheet fantoom: 106764

Note: Alternate identification labels have been supplied with this phantom assembly which provide phantom-handling instructions in a variety of languages. Please select the appropriate label and apply it directly on top of the existing label.

Multi-language phantom label sheet: 106764

HUOMAA: Tämän fantomikokoonpanon mukana on toimitettu erilaisia tunnusmerkkejä, joista löytyvät fantomin käsittelyohjeet eri kielillä. Valitse sopiva merkki ja kiinnitä se suoraan laitemerkin päälle.

Monikielinen fantomin merkkilappu: 106764

REMARQUE : d'autres étiquettes d'identification ont été fournies avec ce fantôme. Elles fournissent des instructions de manipulation du fantôme en plusieurs langues. Sélectionner l'étiquette appropriée et l'appliquer directement par-dessus l'étiquette existante.

Fiche de l'étiquette du fantôme multilingue : 106764

HINWEIS: Mit diesem Phantom wurden Wechseletiketten geliefert. Sie enthalten die Anleitung zur Handhabung des Phantoms in verschiedenen Sprachen. Bitte wählen Sie das passenden Etikett und überkleben Sie damit das vorhandene Etikett.

Mehrsprachiger Phantom-Etikettenbogen: 106764

ΣΗΜΕΙΩΣΗ: Αυτή η διάταξη phantom διατίθεται με εναλλακτικές ετικέτες επισημάνσης, οι οποίες παρέχουν οδηγίες χειρισμού του phantom σε διάφορες γλώσσες. Επιλέξτε την κατάλληλη ετικέτα και τοποθετήστε την απευθείας επάνω από την υπάρχουσα ετικέτα.

Φύλλο ετικετών phantom πολλαπλών γλωσσών: 106764

NOTA: Sono incluse etichette che forniscono le istruzioni di utilizzo in diverse lingue. Selezionare l'etichetta adeguata e applicarla direttamente sopra all'etichetta esistente.

Etichette per fantoccio in lingue diverse: 106764

注意： 代替識別ラベルが、様々な言語で記述されたファントム使用説明書付きのファントムアセンブリと共に提供されます。適切なラベルを選択し、既存のラベルの上に直接貼ってください。

多国言語対応ファントムラベルシート：106764

참고： 이 팬텀 어셈블리와 함께 제공된 대체 ID 레이블에서는 팬텀 취급 지침을 다양한 언어로 전달합니다. 알맞은 레이블을 선택하여 기존 레이블 바로 위에 부착하십시오.

다국어 팬텀 레이블 시트: 106764

MERK: Det medfølger alternative identifikasjonsmerker med denne fantomenheten, som inneholder instruksjoner for fantomhåndtering på flere språk. Velg det passende merket, og påfør det direkte oppå det eksisterende merket.

Flerspråklig fantommerkeark: 106764

UWAGA: Z niniejszym fantomowym zestawem zastępczym dostarczane są alternatywne etykiety identyfikacyjne, zawierające instrukcje w różnych językach. Proszę wybrać stosowną etykietę i umieścić ją na wierzchu istniejącej na zestawie.

Arkusz z etykietami zastępczymi w wielu językach: 106764

NOTA: Juntamente com este fantoma são fornecidos rótulos de identificação alternativos, contendo instruções para a manipulação do fantoma, em vários idiomas. É favor seleccionar o rótulo apropriado e aplicá-lo directamente por cima do rótulo existente.

Folha com os rótulos de fantoma em vários idiomas: 106764

OBSERVAȚIE: etichetele de identificare alternative au fost livrate împreună cu acest ansamblu al fantomei și conțin instrucțiuni de manevrare în diferite limbi. Selectați eticheta corespunzătoare și aplicați-o direct peste cea existentă.

Fișă multilingvistică de etichetare a fantomei: 106764

ПРИМЕЧАНИЕ. С этим фантомным блоком поставляются запасные идентификационные наклейки, которые содержат инструкции по работе с фантомом на разных языках. Выберите соответствующую наклейку и наклейте ее поверх инструкций на блоке.

Таблица многоязычных наклеек для фантома: 106764

Nota: Con este fantoma se suministran etiquetas alternativas de identificación con instrucciones para la manipulación del fantoma en varios idiomas. Por favor seleccione la etiqueta que corresponda y aplíquela directamente sobre la etiqueta existente.

Hoja multilingüe de etiquetas del fantoma: 106764

OBS: De alternativa identifikationsetiketter som medföljer denna fantomenhet tillhandahåller fantomhanteringsinstruktioner på flera olika språk. Välj lämplig etikett och placera den direct ovanpå befintlig etikett.

Ark med etiketter på flera språk: 106764

NOT: Bu fantom tertibatla birlikte çeşitli dillerde fantom kullanma talimatı veren alternatif tanımlama etiketleri sağlanmıştır. Lütfen uygun etiketi seçin ve doğrudan mevcut etiketin üzerine uygulayın.

Çeşitli dillerde fantom etiketi föyü: 106764

pREVISION HISTORY

Revision	Date	Change Description	DCN
01	02-04-06	Initial Issue	N/A
02	03-23-06	GE Part Number was 5145567-2	N/A
03	05-31-06	Addition of Swedish Language Manual Added Swedish Note in Addendum Revised Section 3-3 Invivo Address was Pewaukee Added Revision History	2039
04	11-09-06	Correct EU Representative, Page 2 Corrected GE part number to read 5145567-3 on Japanese version only	2297
05		Added Body Parts Loop Safety Information in Sections 2-6 and 4-3.	2789
06		Added Korean warning note in Addendum	2992
07	02-19-08	Add 9 languages to Operator's Manual package. Add language identifier in top right corner of each translated cover page.	3648

