GE Healthcare

Technical Publications

Direction 5154646-100
Revision 3

CADstream 3.1/4.0 to 4.1 Upgrade Installation/Service Manual

Copyright © 2006 by General Electric Company Inc.
All Rights Reserved

DAMAGE IN TRANSPORTATION

All packages should be closely examined at time of delivery. If damage is apparent, have notation “**damage in shipment**” written on **all** copies of the freight or express bill **before** delivery is accepted or “signed for” by a General Electric representative or a hospital receiving agent. Whether noted or concealed, damage **MUST** be reported to the carrier **immediately** upon discovery, or in any event, within **14** days after receipt, and the contents and containers held for inspection by the carrier. A transportation company will not pay a claim for damage if an inspection is not requested within this **14** day period.

File a report with

- Call 1–800–548–3366 and use option 8.
- Contact your local service coordinator for more information on this process.

Rev. 08/15/2003

WARNING

- THIS SERVICE MANUAL IS AVAILABLE IN ENGLISH ONLY.
- IF A CUSTOMER'S SERVICE PROVIDER REQUIRES A LANGUAGE OTHER THAN ENGLISH, IT IS THE CUSTOMER'S RESPONSIBILITY TO PROVIDE TRANSLATION SERVICES.
- DO NOT ATTEMPT TO SERVICE THE EQUIPMENT UNLESS THIS SERVICE MANUAL HAS BEEN CONSULTED AND IS UNDERSTOOD.
- FAILURE TO HEED THIS WARNING MAY RESULT IN INJURY TO THE SERVICE PROVIDER, OPERATOR OR PATIENT FROM ELECTRIC SHOCK, MECHANICAL OR OTHER HAZARDS.

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REVISION HISTORY

REV	DATE	PRIMARY REASON FOR CHANGE
1.....	December 2, 2005	Initial Version
2.....	February 21, 2006	M3 release
3.....	July10, 2006	M4 release

* This revision/letter corresponds to the indicated document's revision control system.

TABLE OF CONTENTS

1- INTRODUCTION	7
1-1 CADstream 4.1 Advanced Application “Clinical” Summary	7
2- PRELIMINARY REQUIREMENTS	7
2-1 CADstream version.....	7
2-2 Database	7
3- INSTALLATION PROCEDURE	7
3-1 CADstream 3.1/4.0 upgrade procedure	7
3-2 CADstream 3.1/4.0 manual-upgrade procedure.....	9
4- INSTALLATION FINISH TASKS	10
4-1 GIB Update	10
5- SPECIAL CONSIDERATIONS	10
5-1 Preferences	10
6- FRU LIST	10

1- INTRODUCTION

The following manual provides the instruction for the upgrade of the CADstream 3.1 or CADstream 4.0 server to CADstream 4.1. The CADstream 4.1 application is designed to run on the separate CADstream server hw and can be accessed remotely from the Advantage Workstations and MR console. The AW platform should be AW 4.2 and higher. The only MR scanner version that have remote access to the CADstream server is HDx.

1-1 CADstream 4.1 Advanced Application “Clinical” Summary

The CADstream 4.1 application is a post-processing package for Magnetic Resonance Imaging (MRI) breast images. The CADstream 4.1 software that is installed on the AW and OC accesses the CADstream server that is installed separately. Using CADstream 4.1 software, the user can define a process list, process and export CADstream results to DICOM nodes. The processed cases can be analyzed and reviewed using the CADalyst application which is the viewer and allows the users to interact with volumes and AngioMap overlays.

2- PRELIMINARY REQUIREMENTS

2-1 CADstream version

The CADstream 3.1 or 4.0 version should be installed on the CADstream server.

2-2 Database

If the upgrade is done to CADstream 3.1, all studies have to be removed from the CADstream database. They will not be compatible with the 4.1 software. The studies on CADstream 4.0 will be preserved.

3- INSTALLATION PROCEDURE

The automatic upgrade procedure will uninstall the current version of CADstream and will install the CADstream 4.1 version.

3-1 CADstream 3.1/4.0 upgrade procedure

1. Login into the CADstream server
User Name: administrator
Password: adw2.0
2. Insert the upgrade CD into the drive
The upgrade process will be started automatically
3. In Norton Anti-Virus window select: Action – Authorize this script and click OK. See figure 1
4. In Norton Anti-Virus window select: Action – Authorize this script and click OK. See figure 2
5. Click Finish on the CADstream setup window
6. Click Yes on restart window
The system will restart. The upgrade is done.

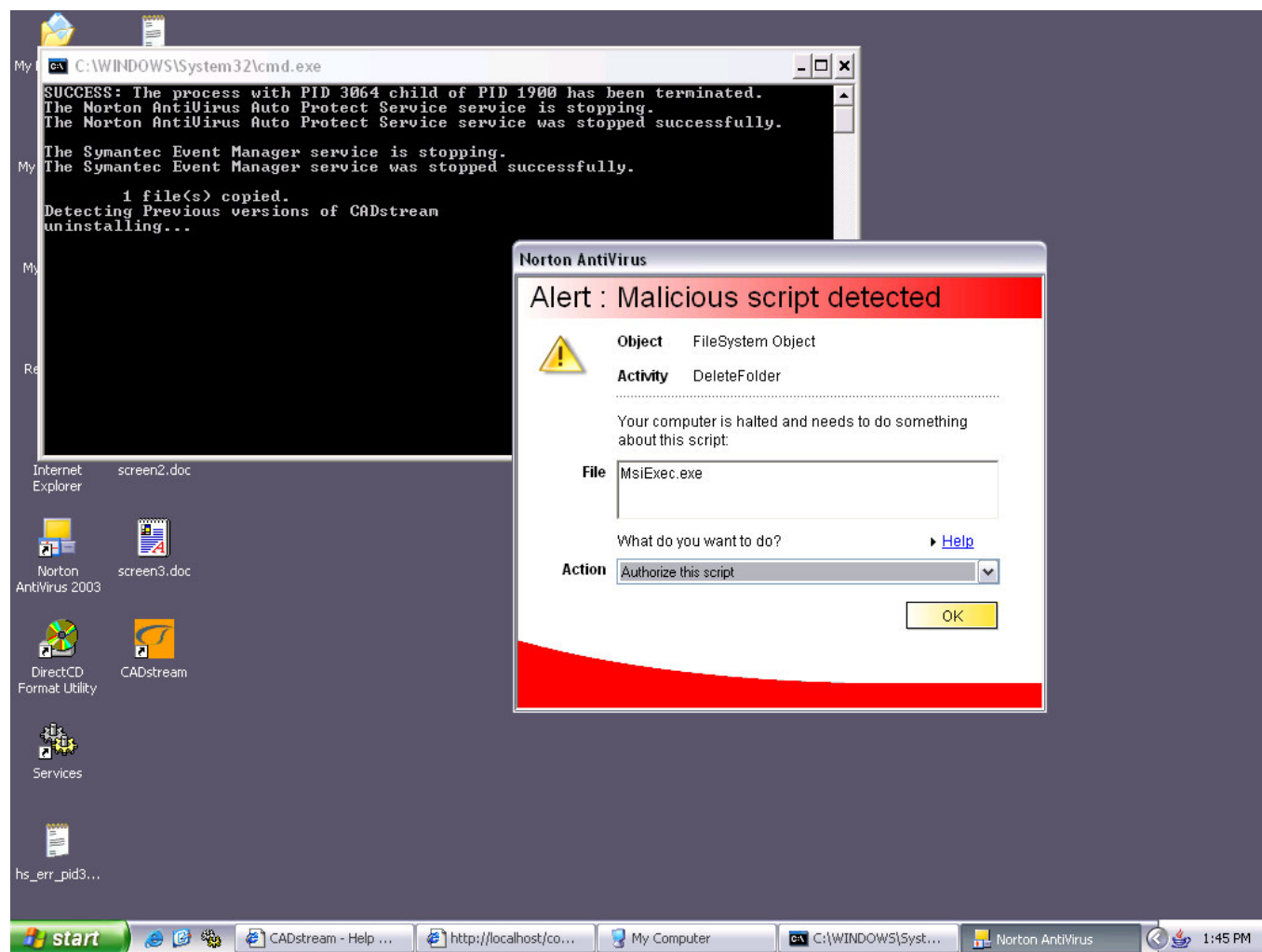
FIGURE 1: FIRST NORTON ANTIVIRUS WINDOW

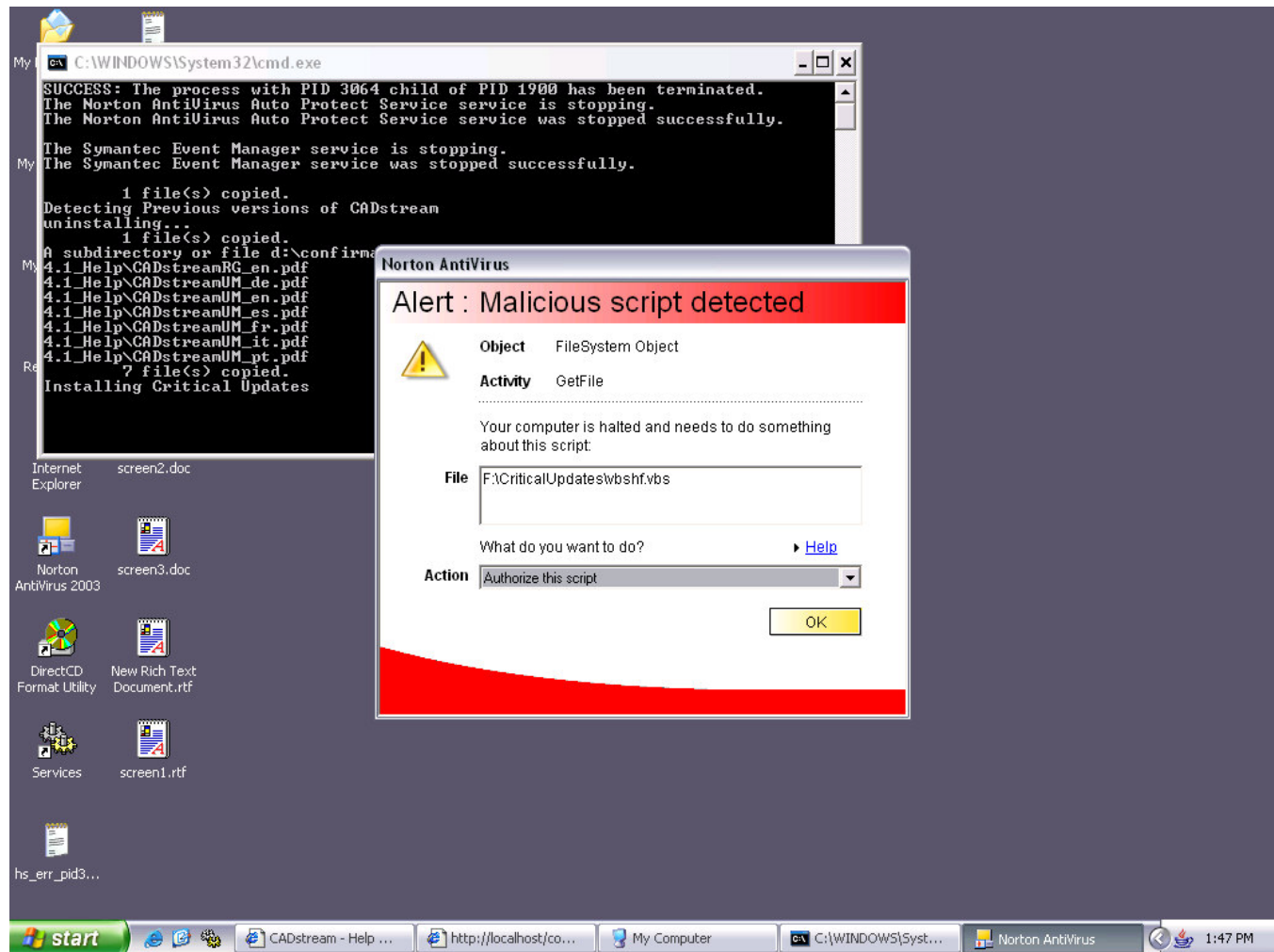
FIGURE 2: SECOND NORTON ANTIVIRUS WINDOW

3-2 CADstream 3.1/4.0 manual-upgrade procedure

1. Should the auto-upgrade fail to complete you will need to manually remove CADstream and re-install the application
2. Login into the CADstream server
User Name: administrator
Password: adw2.0
3. Insert the upgrade CD into the drive, holding down the "shift" key to prevent the auto-upgrade from automatically running
4. Go to the "Start" button, select "Control Panel", and "Add or Remove Programs"
5. Select "CADstream" from the list and chose "Remove". Accept all default prompts during this process. DO NOT restart when prompted.
6. Return to "Add or Remove Programs" and select "Java 2 Runtime Environment" and chose "Remove". Select all default prompts during the uninstall. Restart the computer upon completion.
7. Once the computer has been restarted login to the system as the local administrator
8. Double click the "My Computer" icon and navigate to the F drive (DVD drive)

9. Open the “JRE” folder and double click “JRESetup.exe”. Accept all defaults during the installation.
10. Select the F (DVD) drive again and double click “Setup.exe”. Accept all defaults during the installation. Ensure that the install path is pointing to D:\ (if it is not change the drive letter for install to D:\). Restart the computer upon successful installation.
11. Upon the restart the upgrade is complete. Login to CADstream to verify a successful upgrade.

4- INSTALLATION FINISH TASKS

4-1 GIB Update

The CADstream server and CADstream AW Client application must be entered into GIB (Global Installed Base) database after the initial successful installation in order to fulfill product-tracking requirements.

IF THIS IS NOT THE CASE – PLEASE COMPLETE THIS ASAP AND ENTER THE AW SYSTEM, CADSTREAM SERVER, AS WELL AS CADSTREAM AW CLIENT’S ICD CARD MODEL AND SERIAL NUMBER INTO GIB!

FAILURE TO ENTER THE SYSTEM AND ALL SOFTWARE APPLICATIONS INTO GIB will result in GE Healthcare being in violation of FDA tracking compliance, and this system will not be in the affectivity list for future upgrades.

The OLC AW Support line can direct you as how to get GIB updated correctly if you don’t know. Here is the global GIB link:

GIB: http://egems.gemedicalsystems.com/gib/gib_entry.jsp

5- SPECIAL CONSIDERATIONS

5-1 Preferences

All site settings and preferences will be preserved during the upgrade.

5-2 CADstream server license

New CADstream 4.1 features like SureLoc and Spectroscopy will be automatically enable after the upgrade. No changes are required to the license on the CADstream server.

6- FRU LIST

The replacement parts for this option are as follows:

Part #5160894 – CADstream 4.1 AW Software and Documentation CD Rom

Part#5164701 – CADstream 4.1 (Confirma) Server w/Software

Part#5154402 – CADstream 4.1 (Confirma) Server Software (No Hardware)