


The PHPS (Patient Handling Power Supply) is made up of two AC-DC power supplies, a control board assembly, a servo amplifier for longitudinal motor control, and a CAN Core Communication Board for communicating the PHPS status to the host. The first power supply is the Magnet Room PS. The Magnet Room PS provides switched DC to the PAL and Bore light, and DC power to the SRI, and the PAC. The second power supply is the Table PS. The Table PS provides DC power to the Dock motor and servo amplifier. The PHPS module will reside on the penetration panel in the Excite III base MR system.

The PHPS has a lighted circuit breaker that lights when power is on. The PHPS also has multiple test points for physically checking the power supply voltages. The I/O status can also be checked via the software Diagnostics. From the **MR Service Desktop**, select the **Diagnostics** tab, then open **Peripherals\PHPS Functional Diagnostic**. Select **I/O Status** and then the **RUN** button. The screen will show the status of the Table DC Supply, Bore Light, Bore Vent, PAL (Patient Alignment Light), Table Motor, and Motor Amplifier. The I/O Status will also provide the power supply values for the 8 volt supply, 5 volt supply, and 3.3 volt supply.

**NOTE:** You must turn on the PAL and Bore Light prior to the **I/O Status** test otherwise the status will show as disabled.

The PHPS diagnostic tests run via the CAN communication. If you have communication problems, run the CAN diagnostics which is selectable from the **Diagnostics** tab under **Peripherals**. The PHPS will show up as Node 1.

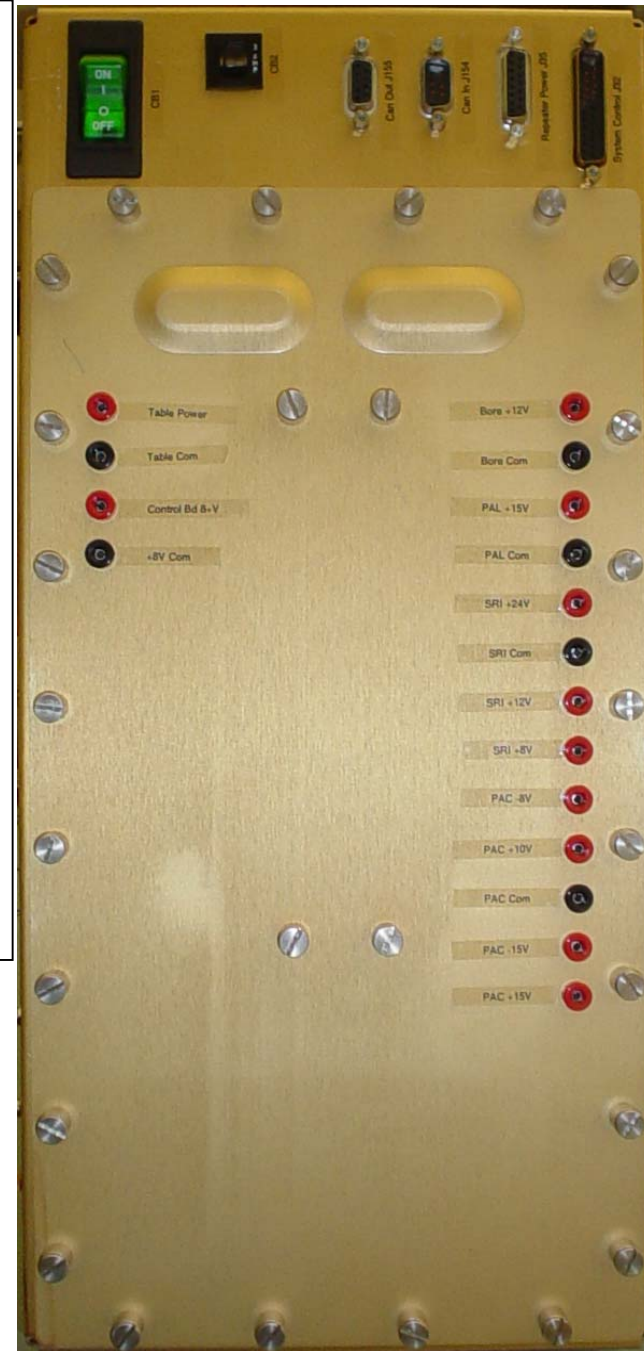

### Available PHPS Diagnostics:

- I/O Status
- Download Amp Config
- Upload Amp Config
- Reset Amp Motor
- Amp Real Time Data