

GE Healthcare

Signa TRM to BRM System Upgrade

Technical Publication



5322720
Rev 4

DAMAGE IN TRANSPORTATION

All packages should be closely examined at time of delivery. If damage is apparent, have notation “**damage in shipment**” written on **all** copies of the freight or express bill **before** delivery is accepted or “signed for” by a General Electric representative or a hospital receiving agent. Whether noted or concealed, damage **MUST** be reported to the carrier **immediately** upon discovery, or in any event, within **14** days after receipt, and the contents and containers held for inspection by the carrier. A transportation company will not pay a claim for damage if an inspection is not requested within this **14** day period.

File a report with

- Call 1–800–548–3366 and use option 8.
- Contact your local service coordinator for more information on this process.

Rev. 08/15/2003

© 2008–2010 General Electric Company
GE Medical Systems, a General Electric Company, going to market as GE Healthcare
3000 N. Grandview Boulevard
Waukesha, Wisconsin 53188
USA

www.gehealthcare.com



ПРЕДУПРЕЖДЕНИЕ

(BG)

- ТОВА УПЪТВАНЕ ЗА РАБОТА Е НАЛИЧНО САМО НА АНГЛИЙСКИ ЕЗИК.
- АКО ДОСТАВЧИКЪТ НА УСЛУГАТА НА КЛИЕНТА ИЗИСКА ЕЗИК, РАЗЛИЧЕН ОТ АНГЛИЙСКИ, ЗАДЪЛЖЕНИЕ НА КЛИЕНТА Е ДА ОСИГУРИ ПРЕВОД.
- НЕ ИЗПОЛЗВАЙТЕ ОБОРУДВАНЕТО ПРЕДИ ДА СТЕ СЕ КОНСУЛТИРАЛИ И РАЗБРАЛИ УПЪТВАНЕТО ЗА РАБОТА.
- НЕСПАЗВАНЕТО НА ТОВА ПРЕДУПРЕЖДЕНИЕ МОЖЕ ДА ДОВЕДЕ ДО НАРАНЯВАНЕ НА ДОСТАВЧИКА НА УСЛУГАТА, ОПЕРАТОРА ИЛИ ПАЦИЕНТ В РЕЗУЛТАТ НА ТОКОВ УДАР ИЛИ МЕХАНИЧНА ИЛИ ДРУГА ОПАСНОСТ.

警告

(ZH-CN)

- 本维修手册仅提供英文版本。
- 如果维修服务提供商需要非英文版本，客户需自行提供翻译服务。
- 未详细阅读和完全理解本维修手册之前，不得进行维修。
- 忽略本警告可能对维修人员，操作员或患者造成触电、机械伤害或其他形式的伤害。

VÝSTRAHA

(CS)

- TENTO PROVOZNÍ NÁVOD EXISTUJE POUZE V ANGLICKÉM JAZYCE.
- V PŘÍPADĚ, ŽE EXTERNÍ SLUŽBA ZÁKAZNÍKŮM POTŘEBUJE NÁVOD V JINÉM JAZYCE, JE ZAJIŠTĚNÍ PŘEKladU DO ODPOVÍDAJÍCÍHO JAZYKA ÚKOLEM ZÁKAZNÍKA.
- NESNAŽTE SE O ÚDRŽBU TOHOTO ZAŘÍZENÍ, ANIŽ BYSTE SI PŘEČETLI TENTO PROVOZNÍ NÁVOD A POCHOPILI JEHO OBSAH.
- V PŘÍPADĚ NEDODRŽOVÁNÍ TÉTO VÝSTRAHY MŮŽE DOJÍT K PORANĚNÍ PRACOVNÍKA PRODEJNÍHO SERVISU, OBSLUŽNÉHO PERSONÁLU NEBO PACIENTŮ Vlivem ELEKTRICKÉHO PROUDU, RESPEKTIVE Vlivem MECHANICKÝCH ČI JINÝCH RIZIK.

ADVARSEL

(DA)

- DENNE SERVICE MANUAL FINDES KUN PÅ ENGELSK.
- HVIS EN KUNDES TEKNIKER HAR BRUG FOR ET ANDET SPROG END ENGELSK, ER DET KUNDENS ANSVAR AT SØRGE FOR OVERSÆTTELSE.
- FORSØG IKKE AT SERVICERE Udstyret medmindre denne service manual har været konsulteret og er forstået.
- MANGLENDE OVERHOLDELSE AF DENNE ADVARSEL KAN MEDFØRE SKADE PÅ GRUND AF ELEKTRISK, MEKANISK ELLER ANDEN FARE FOR TEKNIKEREN, OPERATØREN ELLER PATIENTEN.

WAARSCHUWING

(NL)

- DEZE ONDERHOUDSHANDLEIDING IS ENKEL IN HET ENGELS VERKRIJGBAAR.
- ALS HET ONDERHOUDSPERSONEEL EEN ANDERE TAAL VEREIST, DAN IS DE KLANT VERANTWOORDELIJK VOOR DE VERTALING ERVAN.
- PROBEER DE APPARATUUR NIET TE ONDERHOUDEN VOORDAT DEZE ONDERHOUDSHANDLEIDING WERD GERAADPLEEGD EN BEGREPEN IS.
- INDIEN DEZE WAARSCHUWING NIET WORDT OPGEVOLGD, ZOU HET ONDERHOUDSPERSONEEL, DE OPERATOR OF EEN PATIËNT GEWOND KUNNEN RAKEN ALS GEVOLG VAN EEN ELEKTRISCHE SCHOK, MECHANISCHE OF ANDERE GEVAREN.
- THIS SERVICE MANUAL IS AVAILABLE IN ENGLISH ONLY.

WARNING

(EN)

- IF A CUSTOMER'S SERVICE PROVIDER REQUIRES A LANGUAGE OTHER THAN ENGLISH, IT IS THE CUSTOMER'S RESPONSIBILITY TO PROVIDE TRANSLATION SERVICES.
- DO NOT ATTEMPT TO SERVICE THE EQUIPMENT UNLESS THIS SERVICE MANUAL HAS BEEN CONSULTED AND IS UNDERSTOOD.
- FAILURE TO HEED THIS WARNING MAY RESULT IN INJURY TO THE SERVICE PROVIDER, OPERATOR, OR PATIENT FROM ELECTRIC SHOCK, OR FROM MECHANICAL OR OTHER HAZARDS.

HOIATUS

(ET)

- KÄESOLEV TEENINDUSJUHEND ON SAADAVAL AINULT INGLISE KEELES.
- KUI KLIENDITEENINDUSE OSUTAJA NÕUAB JUHENDIT INGLISE KEELEST ERINEVAS KEELES, VASTUTAB KLIENT TÕLKETEENUSE OSUTAMISE EEST.
- ÄRGE ÜRITAGE SEADMEID TEENINDADA ENNE EELNEVALT KÄESOLEVA TEENINDUSJUHENDIGA TUTVUMIST JA SELLEST ARU SAAMIST.

VAROITUS

(FI)

- KÄESOLEVA HOIATUSE EIRAMINE VÕIB PÕHJUSTADA TEENUSEOSUTAJA, OPERAATORI VÕI PATSIENDI VIGASTAMIST ELEKTRILÖÖGI, MEHAANILISE VÕI MUU OHU TAGAJÄRJEL.
- TÄMÄ HUOLTO-OHJE ON SAATAVILLA VAIN ENGLANNIKSI.
- JOS ASIAKKAAN HUOLTOHENKILÖSTÖ VAATII MUUTA KUIN ENGLANNINKIELISTÄ MATERIAALIA, TARVITTAVAN KÄÄNNÖKSEN HANKKIMINEN ON ASIAKKAAN VASTUULLA.
- ÄLÄ YRITÄ KORJATA LAITTEISTOA ENNEN KUIN OLET VARMASTI LUKENUT JA YMMÄRTÄNYT TÄMÄN HUOLTO-OHJEEN.
- MIKÄLI TÄTÄ VAROITUSTA EI NOUDATETA, SEURAUKSENA VOI OLLA HUOLTOHENKILÖSTÖN, LAITTEISTON KÄYTTÄJÄN TAI POTILAAN VAHINGOITTUMINEN SÄHKÖISKUN, MEKAANISEN VIAN TAI MUUN VAARATILANTEEN VUOKSI.

ATTENTION

(FR)

- CE MANUEL DE SERVICE N'EST DISPONIBLE QU'EN ANGLAIS.
- SI LE TECHNICIEN DU CLIENT A BESOIN DE CE MANUEL DANS UNE AUTRE LANGUE QUE L'ANGLAIS, C'EST AU CLIENT QU'IL INCOMBE DE LE FAIRE TRADUIRE.
- NE PAS TENTER D'INTERVENIR SUR LES EQUIPEMENTS TANT QUE LE MANUEL SERVICE N'A PAS ETE CONSULTE ET COMPRIS
- LE NON-RESPECT DE CET AVERTISSEMENT PEUT ENTRAÎNER CHEZ LE TECHNICIEN, L'OPÉRATEUR OU LE PATIENT DES BLESSURES DUES À DES DANGERS ÉLECTRIQUES, MÉCANIQUES OU AUTRES.

WARNUNG

(DE)

- DIESE SERVICEANLEITUNG EXISTIERT NUR IN ENGLISCHER SPRACHE.
- FALLS EIN FREMDER KUNDENDIENST EINE ANDERE SPRACHE BENÖTIGT, IST ES AUFGABE DES KUNDEN FÜR EINE ENTSPRECHENDE ÜBERSETZUNG ZU SORGEN.
- VERSUCHEN SIE NICHT DIESE ANLAGE ZU WARTEN, OHNE DIESE SERVICEANLEITUNG GELESEN UND VERSTANDEN ZU HABEN.
- WIRD DIESE WARNUNG NICHT BEACHTET, SO KANN ES ZU VERLETZUNGEN DES KUNDENDIENSTTECHNIKERS, DES BEDIENERS ODER DES PATIENTEN DURCH STROMSCHLÄGE, MECHANISCHE ODER SONSTIGE GEFAHREN KOMMEN.

ΠΡΟΕΙΔΟΠΟΙΗΣΗ

(EL)

- ΤΟ ΠΑΡΟΝ ΕΓΧΕΙΡΙΔΙΟ ΣΕΡΒΙΣ ΔΙΑΤΙΘΕΤΑΙ ΣΤΑ ΑΓΓΛΙΚΑ ΜΟΝΟ.
- ΕΑΝ ΤΟ ΑΤΟΜΟ ΠΑΡΟΧΗΣ ΣΕΡΒΙΣ ΕΝΟΣ ΠΕΛΑΤΗ ΑΠΑΙΤΕΙ ΤΟ ΠΑΡΟΝ ΕΓΧΕΙΡΙΔΙΟ ΣΕ ΓΛΩΣΣΑ ΕΚΤΟΣ ΤΩΝ ΑΓΓΛΙΚΩΝ, ΑΠΟΤΕΛΕΙ ΕΥΘΥΝΗ ΤΟΥ ΠΕΛΑΤΗ ΝΑ ΠΑΡΕΧΕΙ ΥΠΗΡΕΣΙΕΣ ΜΕΤΑΦΡΑΣΗΣ.
- ΜΗΝ ΕΠΙΧΕΙΡΗΣΕΤΕ ΤΗΝ ΕΚΤΕΛΕΣΗ ΕΡΓΑΣΙΩΝ ΣΕΡΒΙΣ ΣΤΟΝ ΕΞΟΠΛΙΣΜΟ ΕΚΤΟΣ ΕΑΝ ΕΧΕΤΕ ΣΥΜΒΟΥΛΕΥΤΕΙ ΚΑΙ ΕΧΕΤΕ ΚΑΤΑΝΟΗΣΕΙ ΤΟ ΠΑΡΟΝ ΕΓΧΕΙΡΙΔΙΟ ΣΕΡΒΙΣ.
- ΕΑΝ ΔΕ ΛΑΒΕΤΕ ΥΠΟΨΗ ΤΗΝ ΠΡΟΕΙΔΟΠΟΙΗΣΗ ΑΥΤΗ, ΕΝΔΕΧΕΤΑΙ ΝΑ ΠΡΟΚΛΗΘΕΙ ΤΡΑΥΜΑΤΙΣΜΟΣ ΣΤΟ ΑΤΟΜΟ ΠΑΡΟΧΗΣ ΣΕΡΒΙΣ, ΣΤΟ ΧΕΙΡΙΣΤΗ Ή ΣΤΟΝ ΑΣΘΕΝΗ ΑΠΟ ΗΛΕΚΤΡΟΠΛΗΞΙΑ, ΜΗΧΑΝΙΚΟΥΣ Ή ΑΛΛΟΥΣ ΚΙΝΔΥΝΟΥΣ.

FIGYELMEZTETÉS

(HU)

- EZEN KARBANTARTÁSI KÉZIKÖNYV KIZÁRÓLAG ANGOL NYELVEN ÉRHEŐ EL.
- HA A VEVŐ SZOLGÁLTATÓJA ANGOLTÓL ELTÉRŐ NYELVRE TART IGÉNYT, AKKOR A VEVŐ FELELŐSSÉGE A FORDÍTÁS ELKÉSZÍTTETÉSE.
- NE PRÓBÁLJA ELKEZDENI HASZNÁLNI A BERENDEZÉST, AMÍG A KARBANTARTÁSI KÉZIKÖNYVBEN LEÍRTAKAT NEM ÉRTELMEZTÉK.
- EZEN FIGYELMEZTETÉS FIGYELMEN KÍVÜL HAGYÁSA A SZOLGÁLTATÓ, MŰKÖDTETŐ VAGY A BETEG ÁRAMÚTÉS, MECHANIKAI VAGY EGYÉB VESZÉLYHELYZET MIATTI SÉRÜLÉSÉT EREDMÉNYEZHETI.

AÐVÖRUN

(IS)

- ÞESSI ÞJÓNUSTUHANDBÓK ER EINGÖNGU FÁANLEG Á ENSKU.
- EF AÐ ÞJÓNUSTUVEITANDI VIÐSKIPTAMANNS ÞARFNAST ANNAS TUNGUMÁLS EN ENSKU, ER ÞAÐ SKYLDA VIÐSKIPTAMANNS AÐ SKAFFA TUNGUMÁLAPJÓNUSTU.
- REYNIÐ EKKI AÐ AFGREIÐA TÆKIÐ NEMA AÐ ÞESSI ÞJÓNUSTUHANDBÓK HEFUR VERIÐ SKOÐUÐ OG SKILIN.
- BROT Á SINNA ÞESSARI AÐVÖRUN GETUR LEITT TIL MEIÐSLA Á ÞJÓNUSTUVEITANDA, STJÓRNANDA EÐA SJÚKLINGS FRÁ RAFLOSTI, VÉLRÆNU EÐA ÖÐRUM ÁHÆTTUM.

AVVERTENZA

(IT)

- IL PRESENTE MANUALE DI MANUTENZIONE E DISPONIBILE SOLTANTO IN INGLESE.
- SE UN ADDETTO ALLA MANUTENZIONE ESTERNO ALLA GEMS RICHIEDE IL MANUALE IN UNA LINGUA DIVERSA, IL CLIENTE E TENUTO A PROVVEDERE DIRETTAMENTE ALLA TRADUZIONE.
- SI PROCEDA ALLA MANUTENZIONE DELL'APPARECCHIATURA SOLO DOPO AVER CONSULTATO IL PRESENTE MANUALE ED AVERNE COMPRESO IL CONTENUTO
- IL NON RISPETTO DELLA PRESENTE AVVERTENZA POTREBBE FAR COMPIERE OPERAZIONI DA CUI DERIVINO LESIONI ALL'ADDETTO ALLA MANUTENZIONE, ALL'UTILIZZATORE ED AL PAZIENTE PER FOLGORAZIONE ELETTRICA, PER URTI MECCANICI OD ALTRI

警告

(JA)

- このサービスマニュアルには英語版しかありません。
- サービスを担当される業者が英語以外の言語を要求される場合、翻訳作業はその業者の責任で行うものとさせていただきます。
- このサービスマニュアルを熟読し理解せずに、装置のサービスを行わないでください。
- この警告に従わない場合、サービスを担当される方、操作員あるいは患者さんが、感電や機械的又はその他の危険により負傷する可能性があります。

경고

(KO)

- 본 서비스 지침서는 영어로만 이용하실 수 있습니다.
- 고객의 서비스 제공자가 영어 이외의 언어를 요구할 경우, 번역 서비스를 제공하는 것은 고객의 책임입니다.
- 본 서비스 지침서를 참고했고 이해하지 않는 한은 해당 장비를 수리하려고 시도하지 마십시오.
- 이 경고에 유의하지 않으면 전기 쇼크, 기계상의 혹은 다른 위험으로부터 서비스 제공자, 운영자 혹은 환자에게 위해를 가할 수 있습니다.

BRĪDINĀJUMS

(LV)

- ŠĪ APKALPES ROKASGRĀMATA IR PIEEJAMA TIKAI ANGLŪ VALODĀ.
- JA KLIENTA APKALPES SNIEDZĒJAM NEPIECIEŠAMA INFORMĀCIJA CITĀ VALODĀ, NEVIS ANGLŪ, KLIENTA PIENĀKUMS IR NODROŠINĀT TULKOŠANU.
- NEVEICIET APRĪKOJUMA APKALPI BEZ APKALPES ROKASGRĀMATAS IZLASĪŠANAS UN SAPRAŠANAS.
- ŠĪ BRĪDINĀJUMA NEIEVĒROŠANA VAR RADĪT ELEKTRISKĀS STRĀVAS TRIECIENA, MEHĀNISKU VAI CITU RISKU IZRAISĪTU TRAUMU APKALPES SNIEDZĒJAM, OPERATORAM VAI PACIENTAM.

ISPĒJIMAS

(LT)

- ŠIS EKSPLOATAVIMO VADOVAS YRA PRIEINAMAS TIK ANGLŪ KALBA.
- JEI KLIENTO PASLAUGŪ TIEKĒJAS REIKALAUJA VADOVO KITA KALBA – NE ANGLŪ, NUMATYTI VERTIMO PASLAUGAS YRA KLIENTO ATSAKOMYBĒ.
- NEMĒGINKITE ATLIKTI ĮRANGOS TECHNINĖS PRIEŽIŪROS, NEBENT ATSIŽVELGĖTE Į ŠĮ EKSPLOATAVIMO VADOVĄ IR JĮ SUPRATOTE.
- JEI NEATKREIPSITE DĖMESIO Į ŠĮ PERSPĖJIMĄ, GALIMI SUŽALOJIMAI DĖL ELEKTROS ŠOKO,
- MECHANINIŲ AR KITŲ PAVOJŲ PASLAUGŲ TIEKĖJUI, OPERATORIUI AR PACIENTUI.

ADVARSEL

(NO)

- DENNE SERVICEHÅNDBOKEN FINNES BARE PÅ ENGELSK.
- HVIS KUNDENS SERVICELEVERANDØR TRENGER ET ANNET SPRÅK, ER DET KUNDENS ANSVAR Å SØRGE FOR OVERSETTELSE.
- IKKE FORSØK Å REPARERE UTSTYRET UTEN AT DENNE SERVICEHÅNDBOKEN ER LEST OG FORSTÅTT.
- MANGLENDE HENSYN TIL DENNE ADVARSELEN KAN FØRE TIL AT SERVICELEVERANDØREN, OPERATØREN ELLER PASIENTEN SKADES PÅ GRUNN AV ELEKTRISK STØT, MEKANISKE ELLER ANDRE FARER.

OSTRZEŻENIE

(PL)

- NINIEJSZY PODRĘCZNIK SERWISOWY DOSTĘPNY JEST JEDYNIIE W JĘZYKU ANGIELSKIM.
- JEŚLI DOSTAWCA USŁUG KLIENTA WYMAGA JĘZYKA INNEGO NIŻ ANGIELSKI, ZAPEWNIENIE USŁUGI TŁUMACZENIA JEST OBOWIĄZKIEM KLIENTA.
- NIE PRÓBOWAĆ SERWISOWAĆ WYPOSAŻENIA BEZ ZAPOZNANIA SIĘ I ZROZUMIENIA NINIEJSZEGO PODRĘCZNIKA SERWISOWEGO.
- NIEZASTOSOWANIE SIĘ DO TEGO OSTRZEŻENIA MOŻE SPOWODOWAĆ URAZY DOSTAWCY USŁUG, OPERATORA LUB PACJENTA W WYNIKU PORAŻENIA ELEKTRYCZNEGO, ZAGROŻENIA MECHANICZNEGO BADŹ INNEGO.

ATENÇÃO

(PT)

- ESTE MANUAL DE ASSISTÊNCIA TÉCNICA SÓ SE ENCONTRA DISPONÍVEL EM INGLÊS.
- SE QUALQUER OUTRO SERVIÇO DE ASSISTÊNCIA TÉCNICA, QUE NÃO A GEMS, SOLICITAR ESTES MANUAIS NOUTRO IDIOMA, É DA RESPONSABILIDADE DO CLIENTE FORNECER OS SERVIÇOS DE TRADUÇÃO.
- NÃO TENHA TENTADO REPARAR O EQUIPAMENTO SEM TER CONSULTADO E COMPREENDIDO ESTE MANUAL DE ASSISTÊNCIA TÉCNICA
- O NÃO CUMPRIMENTO DESTES AVISOS PODE POR EM PERIGO A SEGURANÇA DO TÉCNICO, OPERADOR OU PACIENTE DEVIDO A CHOQUES ELÉTRICOS, MECÂNICOS OU OUTROS.

ATENȚIE

(RO)

- ACEST MANUAL DE SERVICE ESTE DISPONIBIL NUMAI ÎN LIMBA ENGLEZĂ.
- DACĂ UN FURNIZOR DE SERVICII PENTRU CLIEȚI NECESITĂ O ALTĂ LIMBĂ DECÂT CEA ENGLEZĂ, ESTE DE DATORIA CLIENTULUI SĂ FURNIZEZE O TRADUCERE.
- NU ÎNCERCAȚI SĂ REPARAȚI ECHIPAMENTUL DECÂT ULTERIOR CONSULTĂRII ȘI ÎNȚELEGERII ACESTUI MANUAL DE SERVICE.
- IGNORAREA ACESTUI AVERTISMENT AR PUTEA DUCE LA RĂNIREA DEPARANATORULUI, OPERATORULUI SAU PACIENTULUI ÎN URMA PERICOLELOR DE ELECTROCUTARE, MECANICE SAU DE ALTĂ NATURĂ.

ОСТОРОЖНО!

(RU)

- ДАННОЕ РУКОВОДСТВО ПО ОБСЛУЖИВАНИЮ ПРЕДЛАГАЕТСЯ ТОЛЬКО НА АНГЛИЙСКОМ ЯЗЫКЕ.
- ЕСЛИ СЕРВИСНОМУ ПЕРСОНАЛУ КЛИЕНТА НЕОБХОДИМО РУКОВОДСТВО НЕ НА АНГЛИЙСКОМ, А НА КАКОМ-ТО ДРУГОМ ЯЗЫКЕ, КЛИЕНТУ СЛЕДУЕТ САМОСТОЯТЕЛЬНО ОБЕСПЕЧИТЬ ПЕРЕВОД.
- ПЕРЕД ОБСЛУЖИВАНИЕМ ОБОРУДОВАНИЯ ОБЯЗАТЕЛЬНО ОБРАТИТЕСЬ К ДАННОМУ РУКОВОДСТВУ И ПОЙМИТЕ ИЗЛОЖЕННЫЕ В НЕМ СВЕДЕНИЯ.
- НЕСООБЛЮЖДЕНИЕ ТРЕБОВАНИЙ ДАННОГО ПРЕДУПРЕЖДЕНИЯ МОЖЕТ ПРИВЕСТИ К ТОМУ, ЧТО СПЕЦИАЛИСТ ПО ОБСЛУЖИВАНИЮ, ОПЕРАТОР ИЛИ ПАЦИЕНТ ПОЛУЧАТ УДАР ЭЛЕКТРИЧЕСКИМ ТОКОМ, МЕХАНИЧЕСКУЮ ТРАВМУ ИЛИ ДРУГОЕ ПОВРЕЖДЕНИЕ.

UPOZORNENIE

(SK)

- TENTO NÁVOD NA OBSLUHU JE K DISPOZÍCII LEN V ANGLIČTINE.
- AK ZÁKAZNÍKOV POSKYTOVATEĽ SLUŽIEB VYŽADUJE INÝ JAZYK AKO ANGLIČTINU, POSKYTNUTIE PREKLADATEĽSKÝCH SLUŽIEB JE ZODPOVEDNOSŤOU ZÁKAZNÍKA.
- NEPOKÚŠAJTE SA O OBSLUHU ZARIADENIA SKÔR, AKO SI NEPREČÍTATE NÁVOD NA OBSLUHU A NEPOROZUMIETE MU.
- ZANEDBANIE TOHTO UPOZORNENIA MÔŽE VYÚSTIŤ DO ZRANENIA POSKYTOVATEĽA SLUŽIEB, OBSLUHUJÚCEJ OSOBY ALEBO PACIENTA ELEKTRICKÝM PRÚDOM, DO MECHANICKÉHO ALEBO INÉHO NEBEZPEČENSTVA.

ATENCION

(ES)

- ESTE MANUAL DE SERVICIO SOLO EXISTE EN INGLES.
- SI ALGUN PROVEEDOR DE SERVICIOS AJENO A GEMS SOLICITA UN IDIOMA QUE NO SEA EL INGLES, ES RESPONSABILIDAD DEL CLIENTE OFRECER UN SERVICIO DE TRADUCCION
- NO SE DEBERA DAR SERVICIO TECNICO AL EQUIPO, SIN HABER CONSULTADO Y COMPRENDIDO ESTE MANUAL DE SERVICIO
- LA NO OBSERVANCIA DEL PRESENTE AVISO PUEDE DAR LUGAR A QUE EL PROVEEDOR DE SERVICIOS, EL OPERADOR O EL PACIENTE SUFRAN LESIONES PROVOCADAS POR CAUSAS ELÉCTRICAS, MECÁNICAS O DE OTRA NATURALEZA.

VARNING

(SV)

- DEN HÄR SERVICEHANDBOKEN FINNS BARA TILLGÄNGLIG PÅ ENGELSKA.
- OM EN KUNDS SERVICETEKNIKER HAR BEHOV AV ETT ANNAT SPRÅK ÄN ENGELSKA ANSVARAR KUNDEN FÖR ATT TILLHANDAHÅLLA ÖVERSÄTTNINGSTJÄNSTER.
- FÖRSÖK INTE UTFÖRA SERVICE PÅ UTRUSTNINGEN OM DU INTE HAR LÄST OCH FÖRSTÅR DEN HÄR SERVICEHANDBOKEN.
- OM DU INTE TAR HÄNSYN TILL DEN HÄR VARNINGEN KAN DET RESULTERA I SKADOR PÅ SERVICETEKNIKERN, OPERATÖREN ELLER PATIENTEN TILL FÖLJD AV ELEKTRISKA STÖTAR, MEKANISKA FAROR ELLER ANDRA FAROR.

DİKKAT

(TR)

- BU SERVİS KILAVUZUNUN SADECE İNGİLİZCESİ MEVCUTTUR.
- EĞER MÜŞTERİ TEKNİSYENİ BU KILAVUZU İNGİLİZCE DIŞINDA BİR BAŞKA LİSANDAN TALEP EDERSE, BUNU TERCÜME ETTİRMEK MÜŞTERİYE DÜŞER.
- SERVİS KILAVUZUNU OKUYUP ANLAMADAN EKİPMANLARA MÜDAHALE ETMEYİNİZ.
- BU UYARIYA UYULMAMASI, ELEKTRİK, MEKANİK VEYA DİĞER TEHLİKELERDEN DOLAYI TEKNİSYEN, OPERATÖR VEYA HASTANIN YARALANMASINA YOL AÇABİLİR.

This page left blank intentionally.

REVISION HISTORY

REV	DATE	PRIMARY REASON FOR CHANGE
-----	------	---------------------------

1	Aug 11, 2008	Initial release for verification.
2	Oct 14, 2008	Initial release for validation.
3	Oct 26, 2008	Updated manual with additional information obtained from validation procedures. Section 1–9: Added step for reconfigure of BRM. Section 2–3: Added step in flowchart for Gradient Chiller Power Down. Section 3: Added note to power down Chiller. Section 4: Added step for removing 8 nylon nuts when removing end bells. Added reference for using IIQ Kit Direction 2318581 for installing bridge support hardware. Revised Section 6: Software Load Procedure.
4	Oct 26, 2010	Sections 1, 2, and 4: Updated Ferrous Material Hazard warning and added 200 gauss warning for movement of component.. Section 4: Updated Blower Box installation handling warnings. Section 5: Added check for proper blower box anchoring. (CAPA 1859293).

LIST OF EFFECTIVE PAGES

<u>PAGE</u>	<u>REV</u>	<u>PAGE</u>	<u>REV</u>	<u>PAGE</u>	<u>REV</u>	<u>PAGE</u>	<u>REV</u>	<u>PAGE</u>	<u>REV</u>
Title Page	4								
Damage in Trans.	–								
Language Policy	‡								
11 to 81	4								

This page left blank intentionally.

TABLE OF CONTENTS

<u>TITLE</u>	<u>PAGE</u>
SECTION 1 – GETTING STARTED	15
1-1 UPGRADE INTRODUCTION	15
1-2 SITE READY CHECK	15
1-3 PRE-DELIVERY FDO CHECKS FOR SIGNA EXCITE	16
1-4 UPGRADE SCHEDULE CONSIDERATIONS	16
1-5 PRODUCT DELIVERY INSTRUCTIONS AND PACKING	16
1-6 DAMAGE IN TRANSPORTATION	16
1-7 DISCARDED MATERIAL RETURN POLICY	16
1-8 PRODUCT LOCATOR	19
1-9 UPGRADE SEQUENCE	19
SECTION 2 – UPGRADE INSTALLATION PROCESS	21
2-1 SITE READINESS	21
2-2 UPGRADE PROCEDURE	21
2-3 SYSTEM MECHANICAL UPGRADE FLOWCHART	22
SECTION 3 – POWER OFF	25
3-1 SYSTEM SHUTDOWN	25
3-2 POWER SHUTDOWN	25
SECTION 4 – TRM TO BRM UPGRADE	27
4-1 DE-INSTALL EXISTING ENCLOSURE COMPONENTS	27
4-2 DE-INSTALL TAC CABINET	45
4-3 DE-INSTALL GRADIENT CABLES	46
4-4 UPGRADE PENETRATION PANEL	48
4-5 INSTALL NEW MAGNET ENCLOSURE COMPONENTS	49
4-6 INSTALL REAR PEDESTAL	62
4-7 INSTALL CHILLER	67
4-8 INSTALL GRADIENT CABLES	69
4-9 FINAL ENCLOSURE INSTALLATION	75
SECTION 5 – UPGRADE COMPLETION	79
5-1 MECHANICAL INSTALLATION COMPLETION	79
5-2 RETURNING SHIPPING EQUIPMENT	79
5-3 MATERIAL DISPOSITION GUIDANCE	79
5-4 FINAL COORDINATION	79
5-5 POWER ON	80
SECTION 6 – LOAD SOFTWARE	81

This page left blank intentionally.

SECTION 1 – GETTING STARTED

1-1 UPGRADE INTRODUCTION

This publication documents the mechanical installation information for upgrading a Signa® 1.5T TwinSpeed with TRM to Signa 1.5T with BRM.


WARNING!

FERROUS MATERIAL HAZARD!

DELIVERY EQUIPMENT AND OTHER TOOLS AND PARTS REQUIRED FOR THIS INSTALLATION CONTAIN FERROUS MATERIAL AND WILL BE STRONGLY ATTRACTED TO MAGNET AND MAY BECOME DANGEROUS PROJECTILES. KEEP ALL FERROUS TOOLS AS FAR FROM THE MAGNET AS POSSIBLE.

NOTICE:

During this hardware upgrade, the magnet may be ramped to full field strength. Individuals working in the magnet room must be MR safety trained.

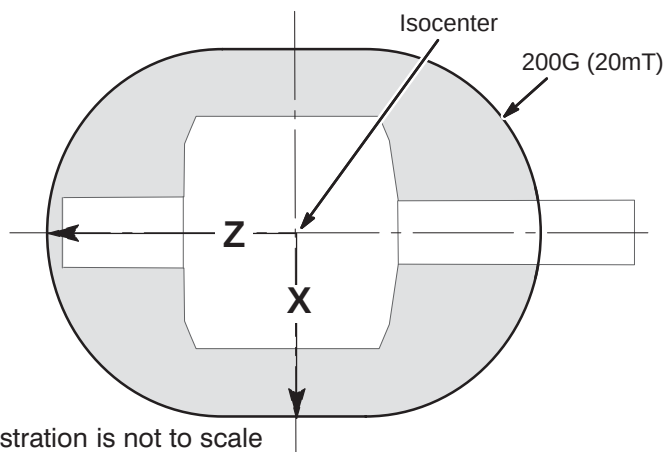

Approximate Distance from isocenter to 200G Line		
	X	Z
1.5T	150 cm (59 in) ¹	205 cm (81 in) ¹
Unshielded Magnet	Measure ²	
Note 1: These distances represent maximum distances for the magnet strength listed. Use gauss meter if needed. Note 2: For unshielded magnets, the magnetic field must be measured using a gauss meter.		

1-2 SITE READY CHECK

Pre-installation work must be completed before equipment is delivered to avoid delays and confusion.

Tools and test equipment to install and calibrate the Signa HDx1.5T system are listed in the *Signa HDx 1.5T Upgrade Pre-Installation* manual, Section 10, TOOLS AND TEST EQUIPMENT.

Service installation tools used in this manual that may not be included in above references lists include:

- 46-301450G1 Fiber Optic Cable Connector Repair kit

1-3 PRE-DELIVERY FDO CHECKS FOR SIGNA EXCITE

Before the shipment is made, Be sure that FDO review has been completed by consultation with the Signa 1.5T Upgrade Specialist. Also verify the Surface coil compatibility issues have also been resolved.

1-4 UPGRADE SCHEDULE CONSIDERATIONS

FMI's, periodic maintenance, installation of prerequisite upgrades and/or options will directly add to the man-hour and system down-time total. The IIS key will need to be re-generated. Be sure that all planned activity is considered when predicting the return of system to customer. Also verify Application Specialist is scheduled at the end of the installation to instruct users. Additional time must also be scheduled to instruct users following completion of upgrade.

1-5 PRODUCT DELIVERY INSTRUCTIONS AND PACKING

The "Shipping Document" lists catalog numbers delivered. **Review and confirm that order is delivered complete.** Check impact on installation schedule if Catalogs and/or packing boxes are missing and/or noted as shipped short. Labels that are attached to the outside of packing boxes summarize box contents. PDI (Product Delivery Instructions) specify box contents, part numbers, and shipping procedure. The PDI is numbered according to the catalog number. Lists of items included with each box are detailed by separate checklists, or a separate sheet that provides a summary of that box contents. Refer to PDI and packing lists for information specific to your shipment.

A set of service and operator manuals is delivered with the Fixed Site Kit. Refer to checklists packed with "Technical Publication" boxes for a list of delivered documents.

1-6 DAMAGE IN TRANSPORTATION

All packages should be closely examined at time of delivery. If damage is apparent, have notation "**damage in shipment**" written on **all** copies of the freight or express bill **before** delivery is accepted or "signed for" by a General Electric representative or a hospital receiving agent. Whether noted or concealed, damage **MUST** be reported to the carrier **immediately** upon discovery, or in any event, within **14** days after receipt, and the contents and containers held for inspection by the carrier. A transportation company will not pay a claim for damage if an inspection is not requested within this **14** day period.

File a report with

- Call 1-800-548-3366 and use option 8.
- Contact your local service coordinator for more information on this process.

1-7 DISCARDED MATERIAL RETURN POLICY

The equipment or material being removed by the upgrade have a potential value for service of the installed base. The disposition of such must be controlled by GE. Disposition of the removed material is per GE Medical Systems Policy and Procedures.

1-7-1 MATERIAL POLICY OVERVIEW

1-7-1-1 Returned Material

Per GE Healthcare Policy 3.9: all requested material is to be returned, as directed, intact and promptly either to GoldSeal or Renewable Resources. Do not remove or swap any parts from the material to be returned.

1-7-1-2 Disposed Material

Regulations governing the disposal of this equipment or material vary from location to location and from country to country. In order to insure compliance with all local, state, or national environment, health, and safety regulations, dispose of all items in the per local GEMS policy and procedures and other applicable published guidance. Contact GEMS Recycling Center for assistance.

1-7-2 MATERIAL TO BE RETURNED

Return all parts, cables, modules, and equipment deinstalled for this upgrade to GoldSeal which includes:

- TRM Body Coil, TRM Gradient Coil with coil cart (if applicable)
- Front and Rear End Bells, and covers
- Comfort Module
- Gradient Cables
- TAC Cabinet
- Coolant Hoses and cables
- TRM Gradient Filter

1-7-3 RETURN KIT

Packaging and/or transportation details can be obtained through the local Installation Specialist. If unavailable, call GoldSeal. The following material may have been provided to facilitate return of material.

- Shipping labels
- Waybill
- Address labels
- Reuse appropriate containers that new hardware and modules were shipped to site for return of above material.
- Packing material

1-7-4 DEINSTALLATION PROCEDURES

Refer to instructions in this manual for deinstallation procedures of site.

1-7-5 REGION SPECIFIC RETURN INFORMATION**1-7-5-1 Sites in Americas**

A GoldSeal return pre-printed BAX waybill and mailing label(s) may be included with the upgrade kit. The upgrade returns kit (p/n 2352804) consists of:

- 2358604 – Waybill
- 2357240 – GoldSeal Return Labels
- 2357241– Upgrade Return Envelope (waybill and return labels)

Gold Seal Address and contact for coordination of transportation arrangements and packaging:

GE Medical Systems – GoldSeal
3114 North Grandview Blvd
Waukesha, WI 53188
Phone: 1-877-251-4468 or 262-513-4147

1-7-5-1 Sites in Americas (Continued)

For non-GoldSeal Returns, utilize GEMS Renewable Resources. For coordination of transportation arrangements and packaging:

Phone: (414)–747–6997 or Dial Comm 8*579–6997

Fax: (414)–747–6855

Field ASPEN: 1–800–525–1516, Box #26041

1-7-5-2 Sites in Asia

At this time, system components are not returned for systems in Asia, but whole systems must be returned from complete system (forklift) upgrades.

Return procedure for GEMS–Asia countries:

Please contact your Regional Service Headquarters to get instructions.

1-7-5-3 Sites in Europe

The GoldSeal Returns Coordinator contacts each site individually and reviews the material to be returned and the logistics for the return. A list of material and a printable shipping label is sent to the site's primary FE, service providers, and forwarders.

GoldSeal Europe return information:

GEMS–E – GOLD SEAL OPERATION & RECYCLING CENTER

GEODIS LOGISTICS – CD 26 – Pièce de la Remise – Bâtiment EVL

Quai 45 à 49 – 91090 LISSES – FRANC

Tel: (+33) 1–60–91–57–35 or (+33) 1 64 97 55 41

Fax: (+33) 1–60–86–12–31

For non-Gold Seal returns in Europe:

European Renewable Resources Center

GEODIS LOGISTICS – CD 26 – Pièce de la Remise – Bâtiment EVL2

Quai 45 à 49 – 91090 LISSES – FRANCE

TEL = (+33)–1–64–97–55–41 FAX = (+33)–1–60–86–12–31

1-7-6 PRODUCT LOCATOR INFORMATION

Process product locator information for all serialized components removed from the system according to policy. Refer to Section 1–8.

1-8 PRODUCT LOCATOR

The Global Install Base Database (also known as Product Locator System) tracks shipment, trans-shipment, and field location of the serialized models. There are now two methods for submitting Product Locator information.

1. At this time, for **U.S. ONLY**, the preferred method of submitting information is the *FE Site Verification Web Site* at: <http://gein2.med.ge.com/gib>

The FE Site Verification consists of three components that are available on the web from the main menu. They are:

- Install/deinstall product locator model and serial numbers
- Add/modify ship to address information
- Update CARES FE data for primary/secondary FEs

To obtain a copy of the FE training tutorial for using this website, a downloadable copy is available at the following: Product Locator Support Central Page – <http://supportcentral.ge.com/15563>

2. One “Shipping Card” is filled out and submitted when shipped (extra cards are supplied for trans-shipments between storage and distribution points), and the “Installation Card” and extra shipment cards are attached.

Verify that serial and model number on each rating plate matches installation card numbers before removing installation card. Note that there may be one or more shipment cards and bar code labels with the installation card. These shipment cards are used to trace the transfer of serialized units between various inventory storage and distribution points until the product reaches its final installation destination. Process just the installation card and discard any extra shipment cards and labels.

1-9 UPGRADE SEQUENCE

Follow the system Upgrade Flow chart in SECTION 2 for best installation efficiency.

Note

All on-site construction must be completed before equipment is delivered and installation starts. Attempting to install the system while construction is being completed will impact installation efficiency and further delay site completion. Making sure that all preinstallation and construction work is completed before equipment is delivered will usually result in an earlier turnover date.

Note that many procedures may be performed in parallel and may be performed in any order according to the specific situation of each site. However all parallel tasks are to be completed prior to the next task.

The general sequence of events is:

- Verify that all required upgrades have arrived complete and site is ready to shut down.
- Shut site down.
- Remove obsolete equipment and cables.

Aside: Removing and installing cables is more efficient and accurate when the equipment room is un-cluttered. It is strongly recommended that this be done after old equipment is removed and before new equipment is positioned. This house cleaning will make the installation tasks much easier.

- Route new cables
- Install new equipment and modify existing equipment.
- All discarded material (i.e. removed/discarded hardware, packing material) must be returned in accordance with GE recycling policy
- Power up
- Load software.
- Perform Functional Checks, System calibrations, and System performance checks
- Reconfigure for BRM
- Replace covers, complete final tasks, and return site to customer.

This page left blank intentionally.

SECTION 2 – UPGRADE INSTALLATION PROCESS

2-1 SITE READINESS

All on-site construction must be completed before equipment is delivered and installation starts. Attempting to upgrade the system while construction is being completed will impact installation efficiency and further delay site completion. Making sure that all preinstallation and construction work is completed before equipment is delivered will usually result in an earlier turnover date.



FERROUS MATERIAL HAZARD!

DELIVERY EQUIPMENT AND OTHER TOOLS AND PARTS REQUIRED FOR THIS UPGRADE CONTAIN FERROUS MATERIAL AND WILL BE STRONGLY ATTRACTED TO MAGNET AND MAY BECOME DANGEROUS PROJECTILES. KEEP ALL FERROUS TOOLS AS FAR FROM THE MAGNET AS POSSIBLE.

2-2 UPGRADE PROCEDURE

The following upgrade flow charts contain information for upgrading to Signa TRM system to a BRM system.

The upgrade flow chart should be followed for an orderly and efficient installation. Note that procedures may be performed in parallel and may be performed in any order according to the specific situation of each site.


This upgrade flowchart has been developed assuming that all Signa HDx equipment has been delivered together. **Make sure that every part required is available before starting.**

2-3 SYSTEM MECHANICAL UPGRADE FLOWCHART

Procedures to be Completed by Field Engineer

MARK OFF
PERFORMED
ITEMS

IMPORTANT : PER ISO
9001 REQUIREMENTS

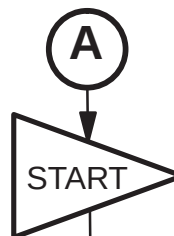

TRM TO BRM UPGRADE CHECKLIST

☐ RE-GENERATE IIS KEY.
This must be done before upgrade starts to
make sure that sytem meets specifications.

☐ PERFORM "saveINFO" PROCEDURE

☐ PERFORM Save Protocol PROCEDURE

Procedures to be Completed by FE and/or Mechanical Installers

Section 1
☐ Section 1-5
PRODUCT DELIVERY INSTRUCTIONS

Section 1
☐ Section 1-6
DAMAGE IN TRANSPORTATION

Section 1
☐ Section 1-8
PRODUCT LOCATOR


Section 3
☐ Section 3-1
SYSTEM POWER SHUTDOWN
☐ LOCKOUT & TAGOUT MAIN DISCONNECT

GRADIENT CHILLER POWER DOWN



2-3 SYSTEM MECHANICAL UPGRADE FLOWCHART (Continued)

De-Installation Procedures to be Completed by FE and/or Mechanical Installers

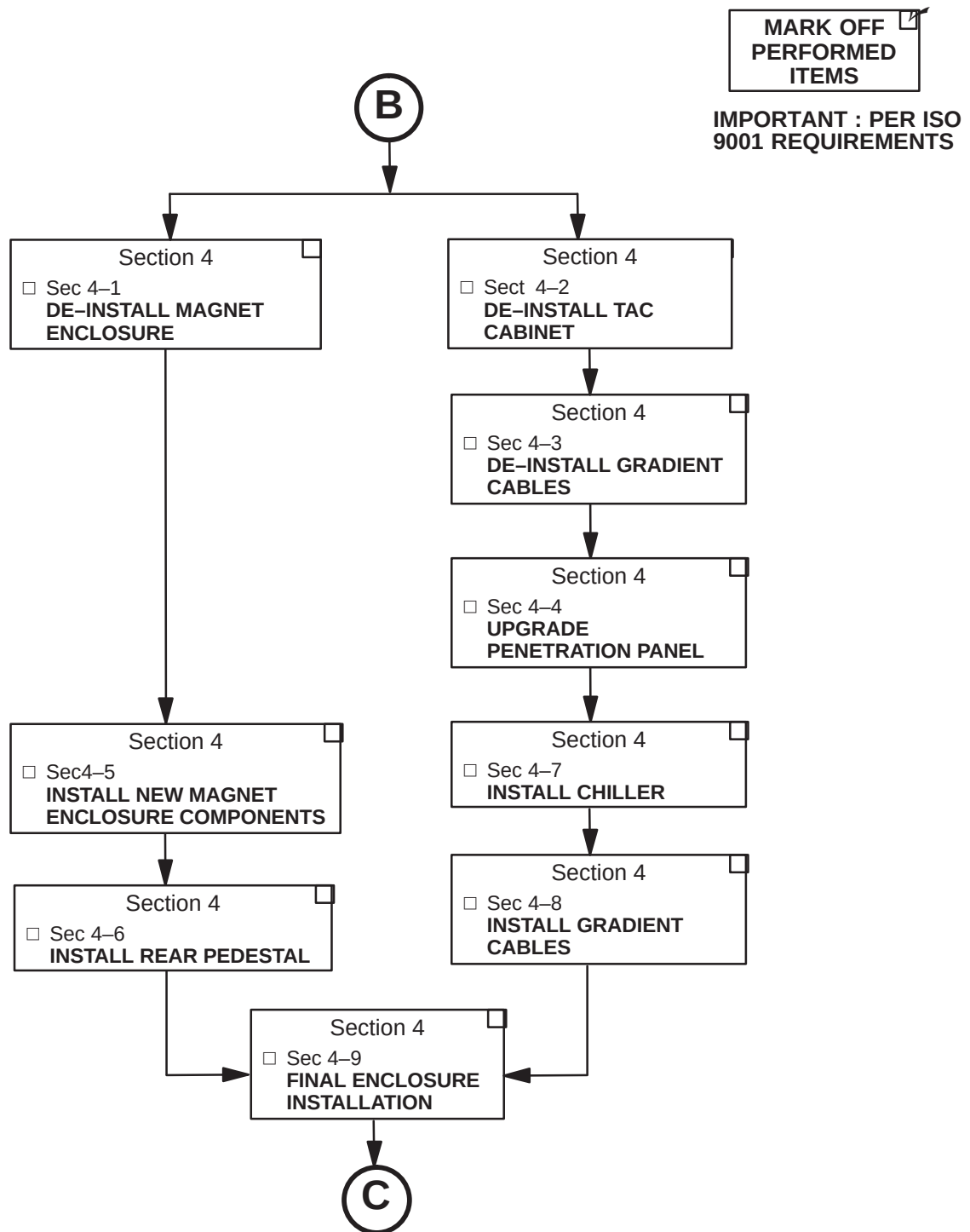

2-3 SYSTEM MECHANICAL UPGRADE FLOWCHART (Continued)

Procedures to be Completed by FE and/or Mechanical Installers

MARK OFF
PERFORMED
ITEMS

IMPORTANT : PER ISO
9001 REQUIREMENTS

Ⓒ

Section 5

- ☐ Section 5-1
**MECHANICAL INSTALLATION
COMPLETION CHECKLIST**
- ☐ Section 5-2
RETURNING SHIPPING EQUIPMENT
- ☐ Section 5-3
MATERIAL DISPOSITION GUIDANCE
- ☐ Section 5-4
FINAL COORDINATION

**MECHANICAL INSTALLATION
COMPLETE**

Procedures to be Completed by FE and Applications

Section 5

- ☐ Section 5-5
POWER ON

MR Service Methods CD-ROM
Direction 5155500

- ☐ Install/Cals>System>Global Operator Console>
LCD Monitor Screen Adjustments

IMPORTANT:
Software MUST BE reloaded.
DO NOT perform Restore Info
as data is for TRM system.
Perform Restore GI only.

Section 6

- ☐ **SOFTWARE LOAD PROCEDURE**
 - MR Service Methods CD-ROM
Software>System Host Software>
Host Operating System Software Load
 - MR Service Methods CD-ROM
Installation>System>
Install and Calibration Wizard

TURN OVER TO CUSTOMER

SECTION 3 – POWER OFF

3-1 SYSTEM SHUTDOWN




FATAL ELECTRIC SHOCK HAZARD!! LETHAL VOLTAGES ARE PRESENT WITHIN THE PDU EVEN IF ALL PDU BREAKERS ARE OFF. TO PREVENT POSSIBLE FATAL ELECTRIC SHOCK VERIFY THAT FACILITY DISCONNECT IS OFF, LOCKED OUT, AND TAGGED BEFORE PROCEEDING. SERIOUS INJURY OR DEATH BY ELECTROCUTION MAY OTHERWISE RESULT.

1. Verify that System performance baseline measurements have been obtained, especially:
 - SPT Full Test Mode
 - White Pixel Test
 - LV Shim
 - IIS
2. Make sure that site customer has archived all patient information. Execute "Save Info" and "Save GI".
3. Perform "Save Protocol" procedure.
4. Record center frequency.
5. Ensure that Patient Transport is undocked with cradle in place and Head Coil carriage is at rear of travel.
6. Perform normal software shutdown procedures.

3-2 POWER SHUTDOWN

1. Press the POWER OFF button on the PDU front control panel.
2. Notify field service and other installation personnel that are working at the site that the PDU main disconnect is to be shut off, locked out, and tagged.
3. Locate and shut off main disconnect supplying power to the PDU.
4. Perform lock out and tag out procedures at main disconnect to prevent inadvertent power restoration.
5. Verify absence of power in the PDU by checking for 0 volts on incoming power wiring

IMPORTANT! – Power down Chiller. The power for the Chiller is supplied by site.

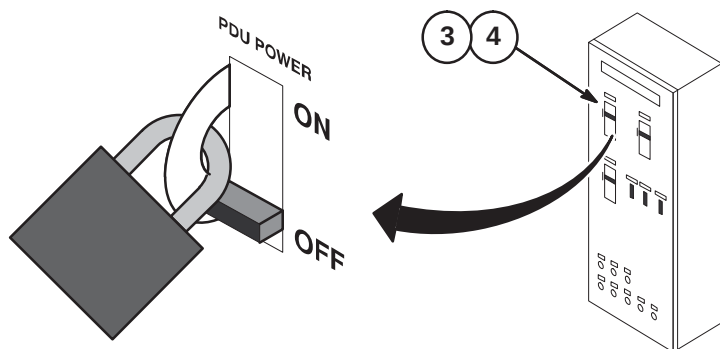


This page left blank intentionally.

SECTION 4 – TRM TO BRM UPGRADE


4-1 DE-INSTALL EXISTING ENCLOSURE COMPONENTS

Before system is powered down, move LPCA to rear of bridge.

4-1-1 Remove Rear Pedestal Cover Plates

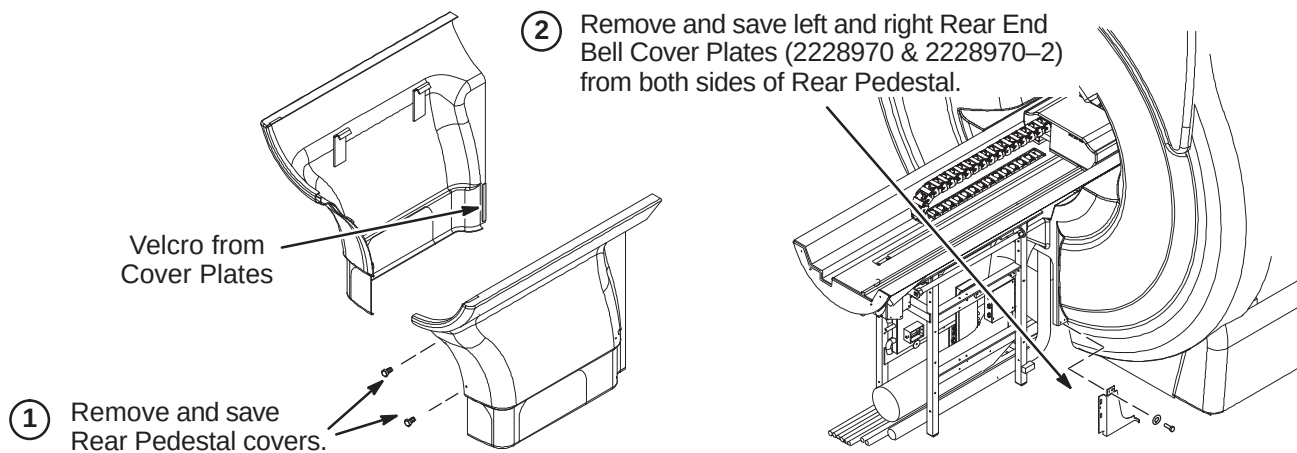


4-1-2 Remove Enclosure Top Covers

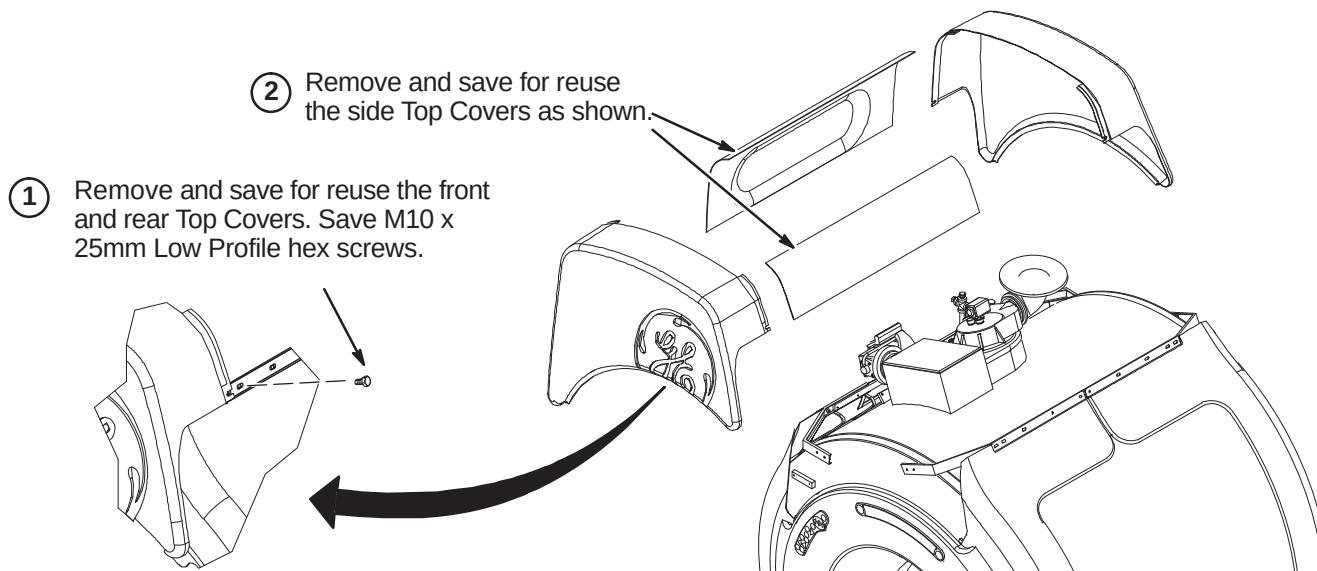



4-1-3 Remove Enclosure Base Trim Covers

Remove and save remaining side Base Trim Covers as shown.

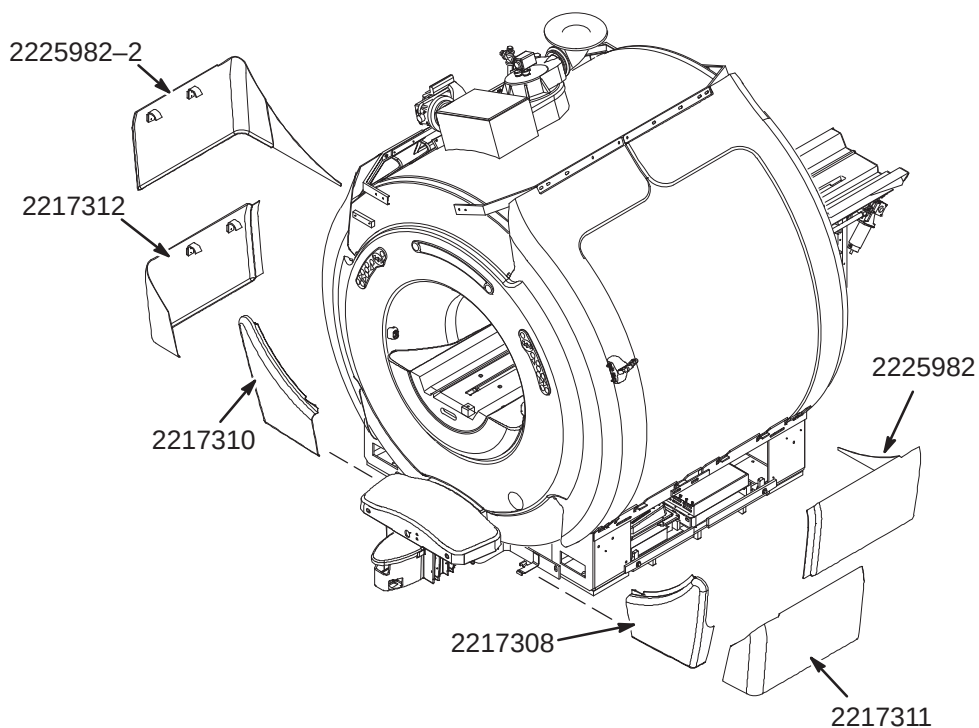


4-1-4 Remove Bridge Cover and Drive Belt

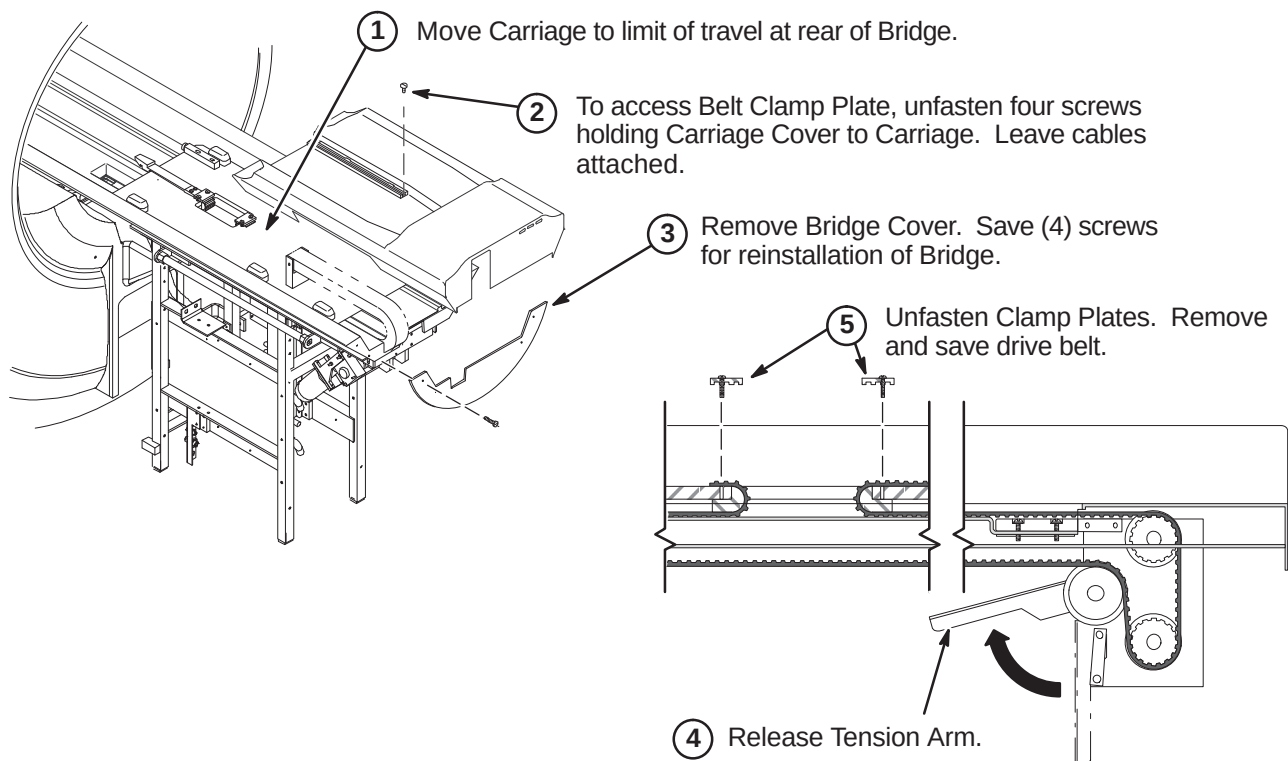

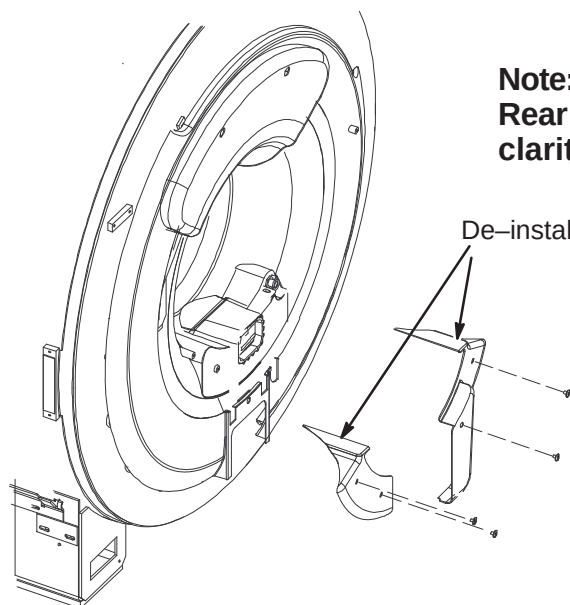
4-1-5 Remove Rear Bridge Covers

Note:
Rear Pedestal and Bridge are not shown for clarity.

De-install and discard left and right Rear Bridge covers and fasteners.

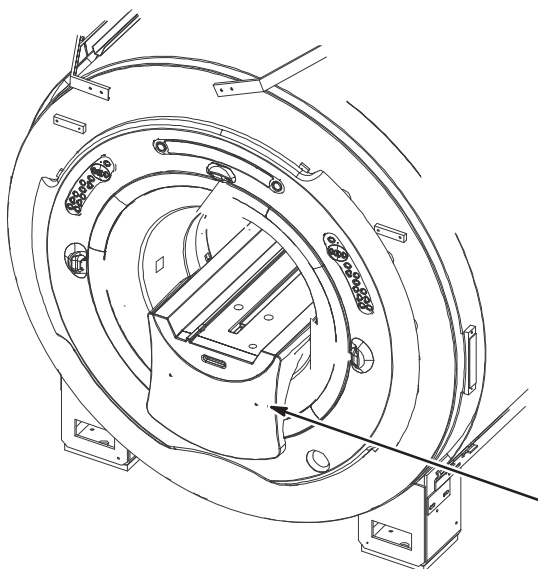
4-1-6 Remove Front Bridge Cover

De-install and discard the Front Bridge Cover and fasteners.

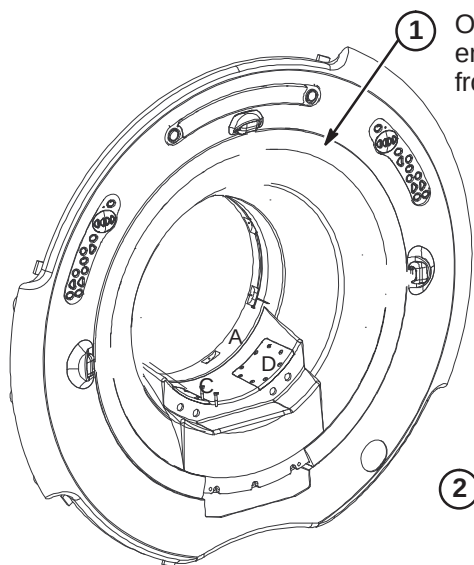
4-1-7 Remove Front and Rear Blue Ring Sections

① On the front of the magnet enclosure, de-install and discard front blue Cosmetic Ring.

② On the rear of the magnet enclosure, de-install and discard two Cosmetic Ring Sections.

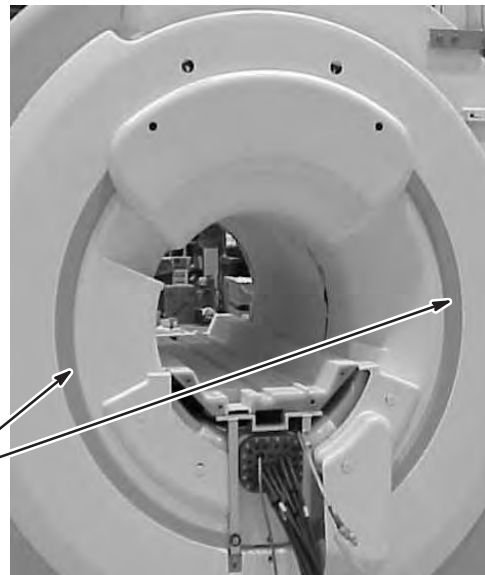**4-1-8 Disconnect Fiber Optic Bundle from Bore Lamp**


Note: Bore Lamp should not have been removed during this upgrade.

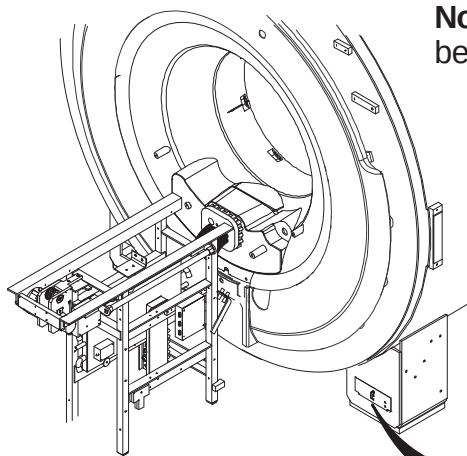


① Bore light Assembly (2222547) should be in place at left rear magnet leg. Carefully remove assembly from opening in leg. The Bore Light assembly will be reused.

③ Remove bushing with inserted fiber optic bundle from clip on Bore Light Assembly.


④ Remove fiber bundle (part of 2222549-4) from Bushing. Save Bushing for re-use. Fiber Optic bundle will be removed with attached light panels

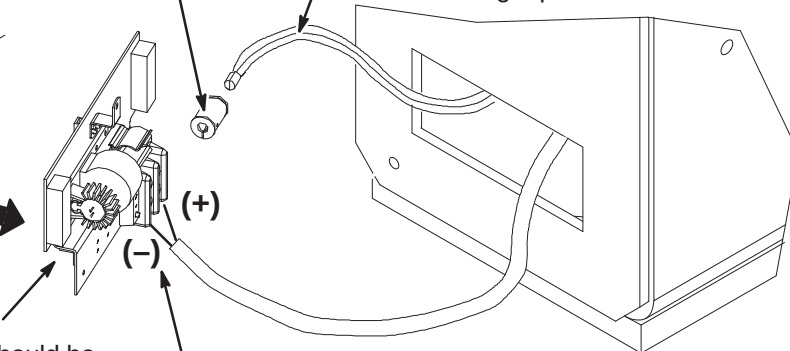


② Do not disconnect the (+) and (-) wires of Run 300 on the terminal strip. Leave connected to terminal strip.

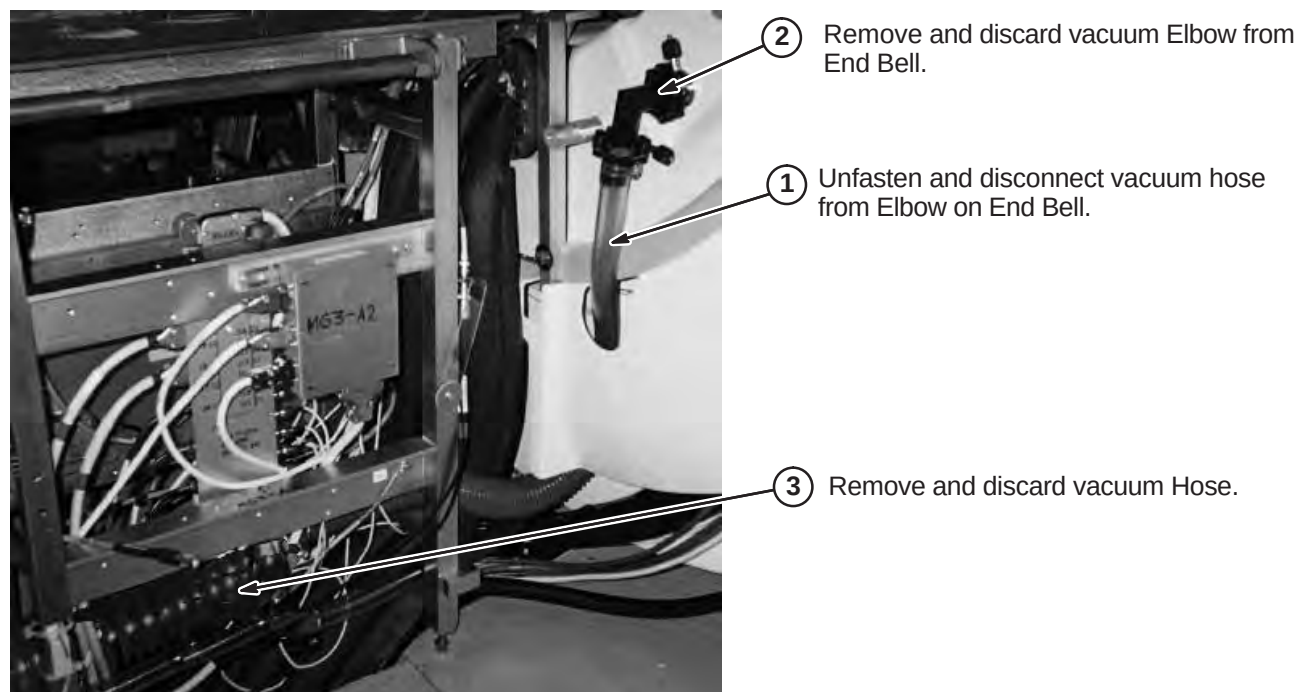4-1-9 Remove Vacuum Hose and Elbow**4-1-10 Disconnect Cables from Magnet to Rear Pedestal**

Note: Cables are disconnected and saved for reconnection to upgraded system.

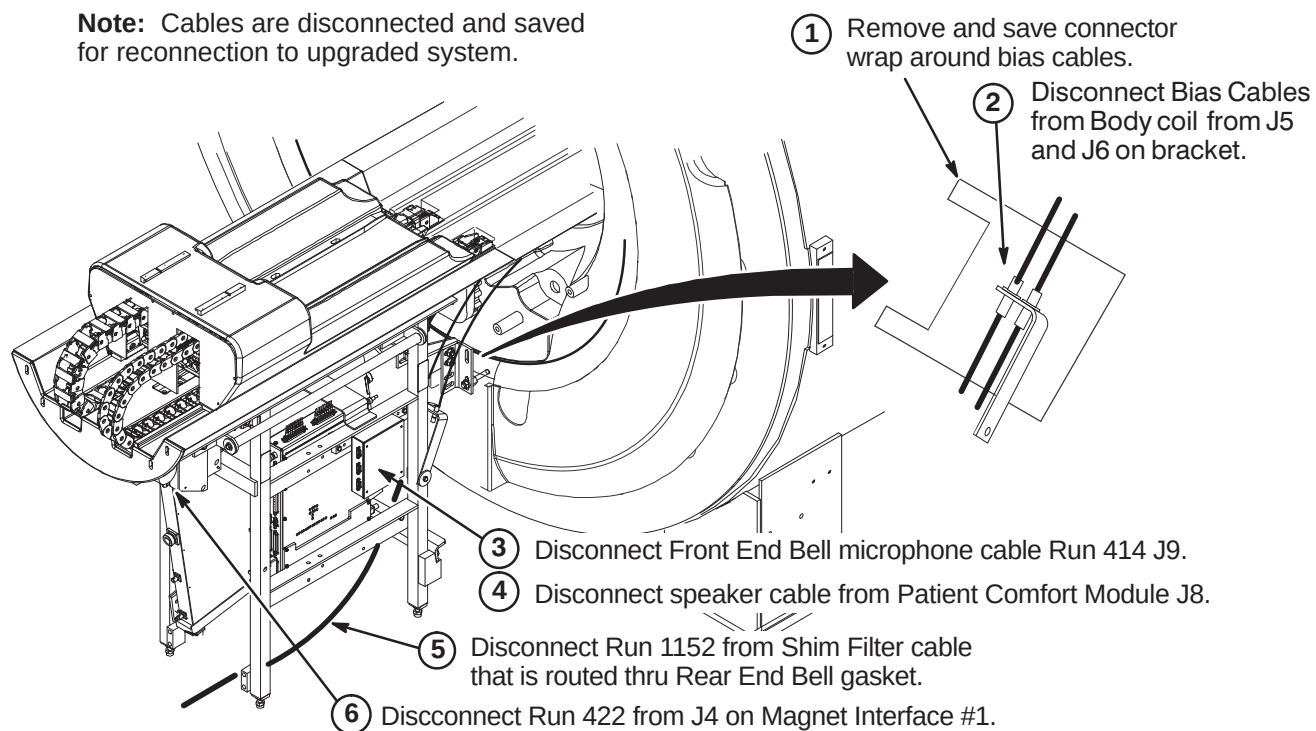

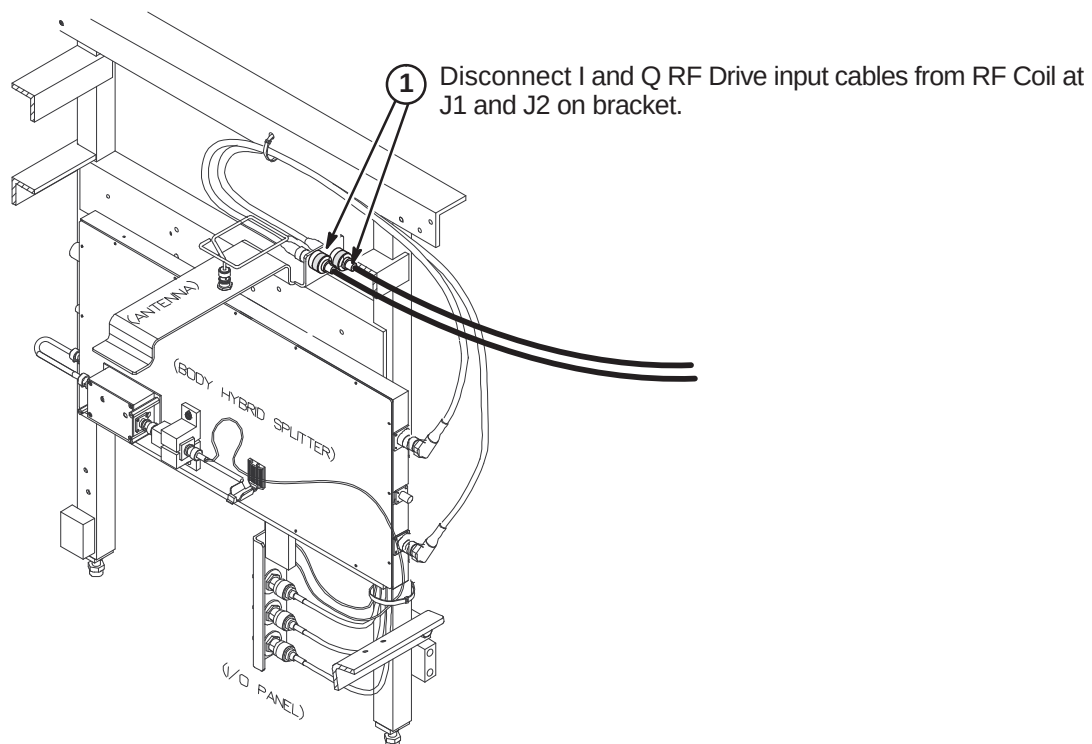
4-1-11 Disconnect I and Q Cables

4-1-12 200 Gauss Zone Limit for Movement of Blower Box


Refer to the illustration and table below for applicable 200 Gauss zone for magnet at site.



THE BLOWER BOX IS HIGHLY FERROUS! IF MAGNET IS RAMPED, ANY MOVEMENT OF THE BLOWER BOX REQUIRES TWO PEOPLE AND MUST BE OUTSIDE THE 200 GAUSS ZONE.

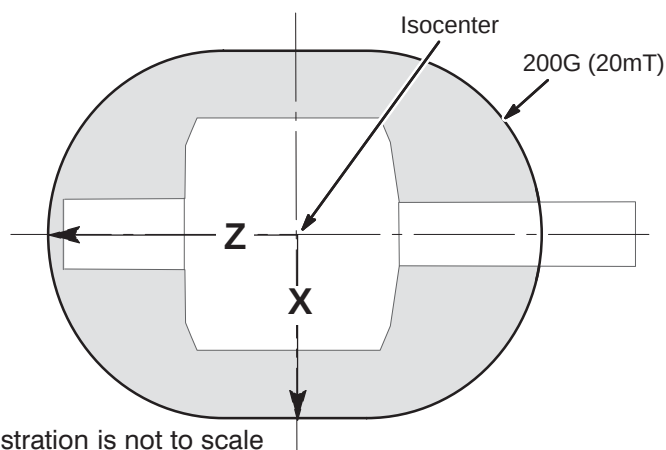


Approximate Distance from isocenter to 200G Line		
	X	Z
1.5T	150 cm (59 in) ¹	205 cm (81 in) ¹
Unshielded Magnet	Measure ²	
Note 1: These distances represent maximum distances for the magnet strength listed. Use gauss meter if needed. Note 2: For unshielded magnets, the magnetic field must be measured using a gauss meter.		


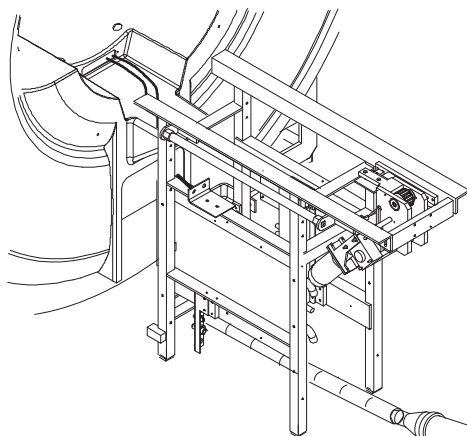
4-1-13 Disconnect and Remove Blower Box



WARNING!

BLOWER BOX CONTAINS FERROUS MATERIAL. REFER TO SECTION 4-1-12 FOR MOVEMENT OF BLOWER BOX.

- ① Disconnect and save 4 inch hose from blower box.
- ② Unfasten blower box mounting hardware and remove. Return for recycling.
- ③ Retain 4 inch hose for reuse in upgraded system.

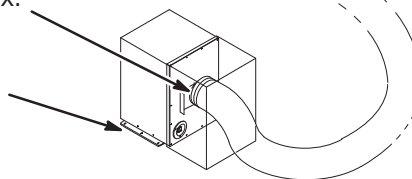

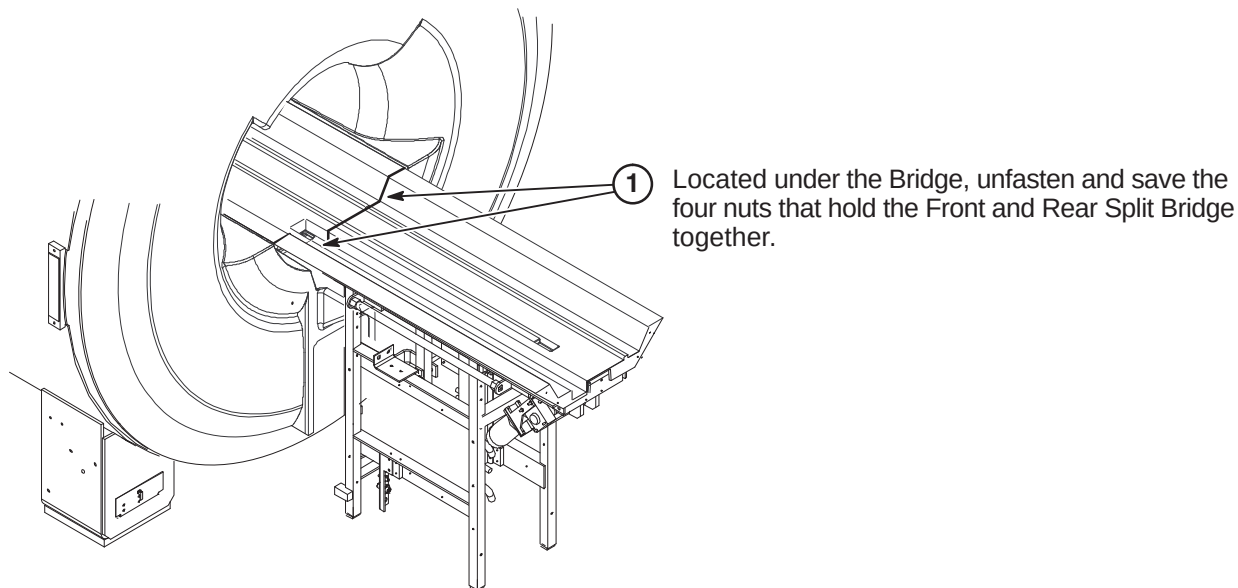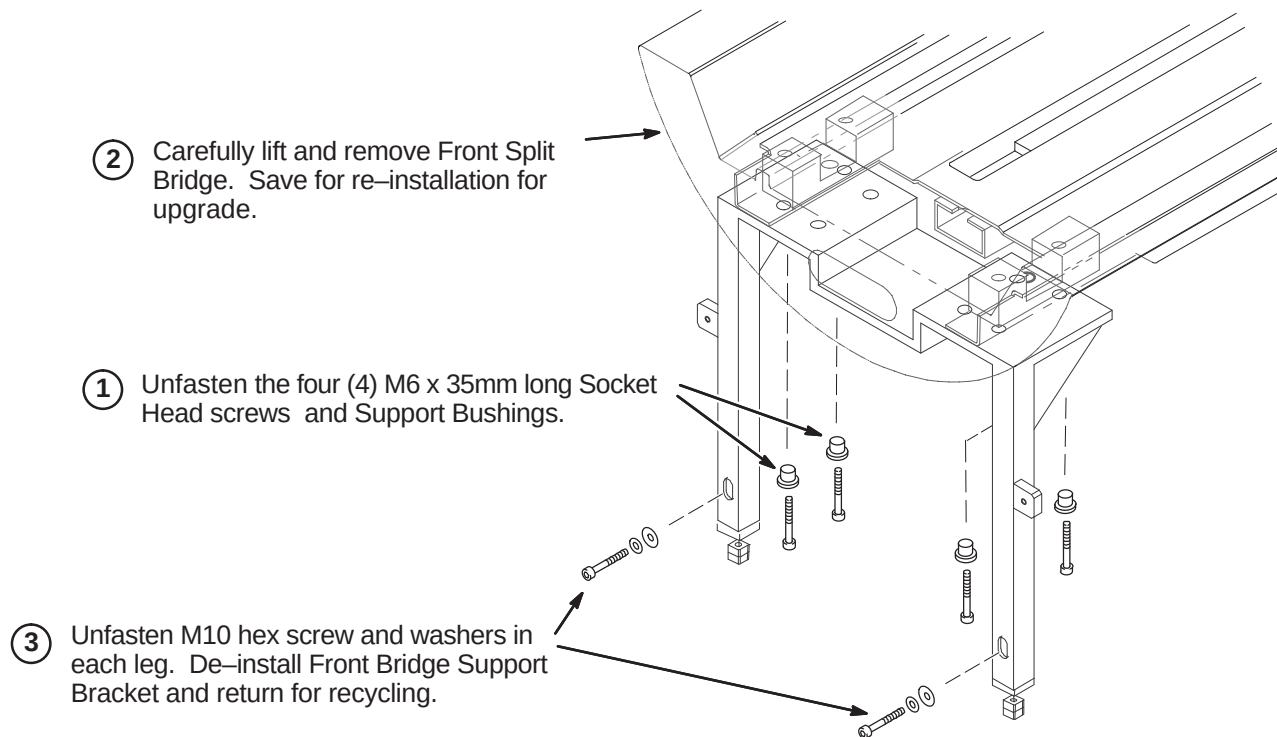4-1-14 Disconnect Front and Rear Split Bridge**4-1-15 Remove Front Split Bridge and Bridge Support**

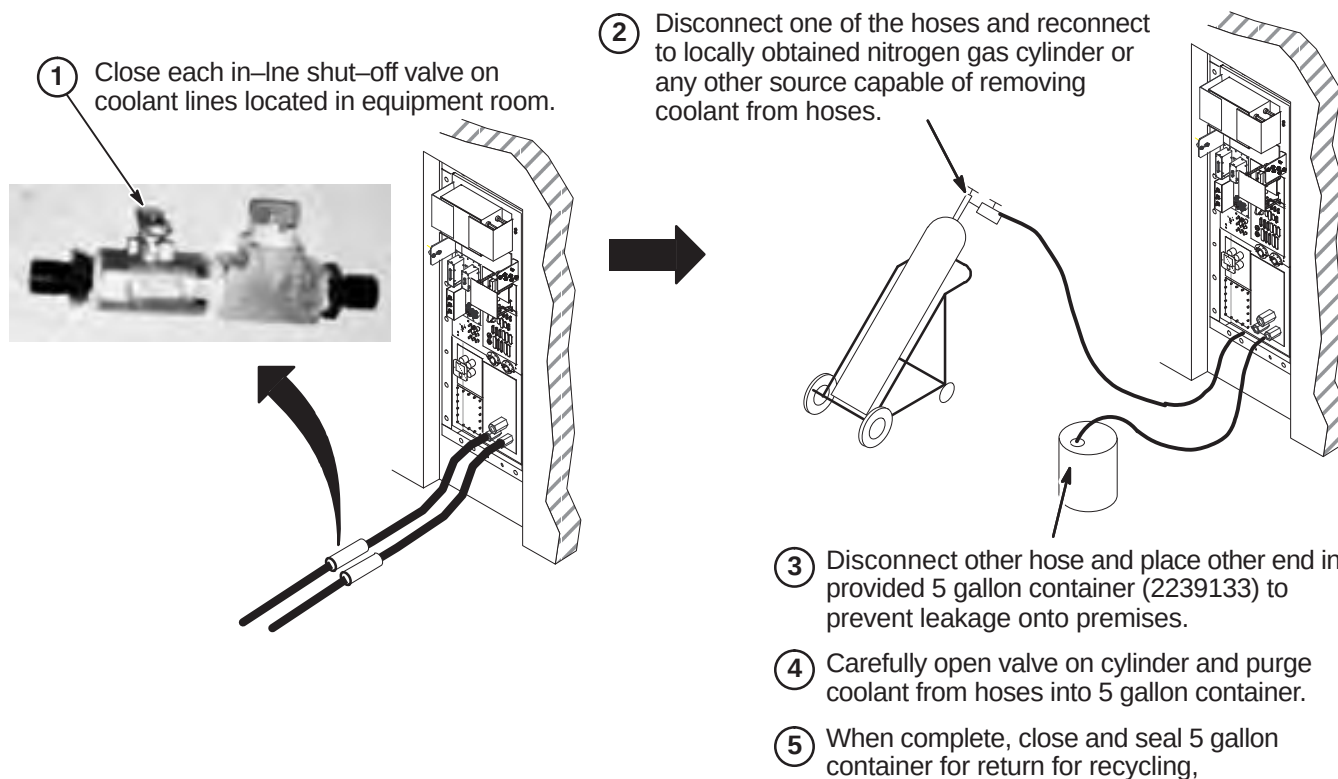
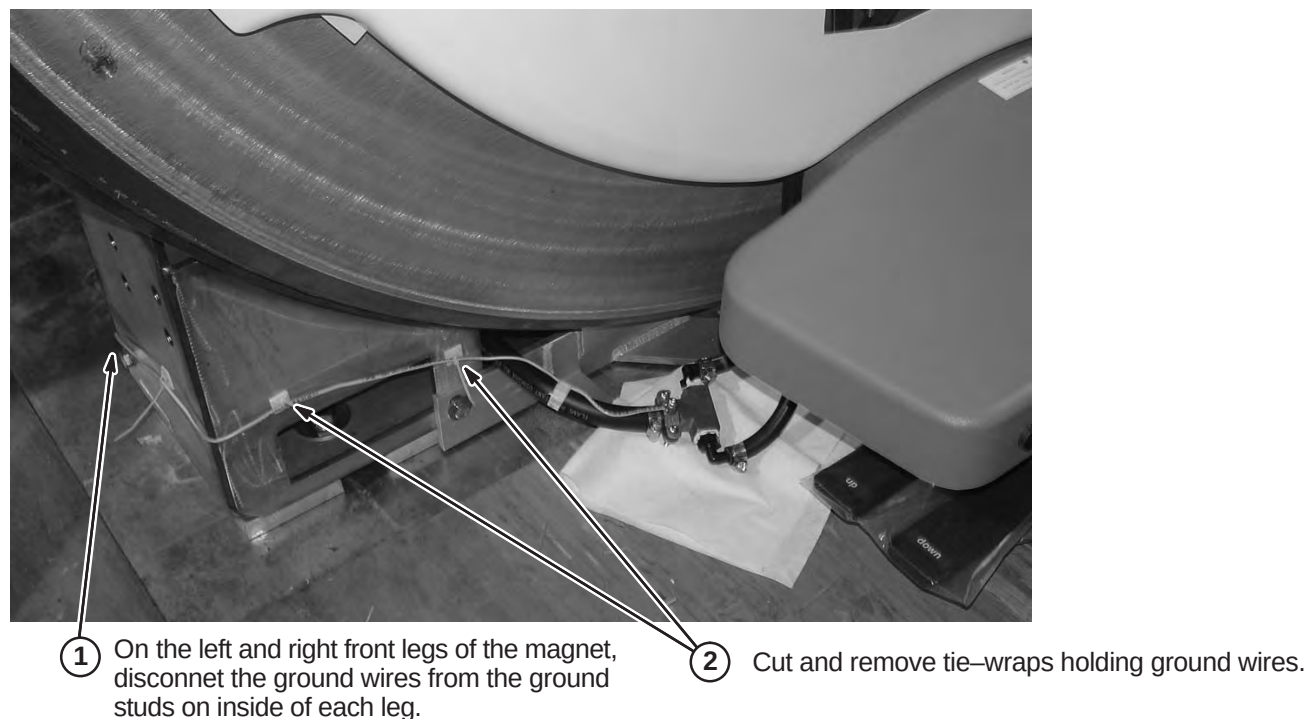
4-1-16 Disconnect Coolant Hoses for Draining**4-1-17 Disconnect Coolant Hose Ground Wires**

4-1-18 Remove Coolant Hoses

① Disconnect the blue and black hoses from hose barbs on end bell.

② Disconnect the blue and black hoses from the 1/2" Union Tees.

③ Remove all coolant hoses and return for recycling.

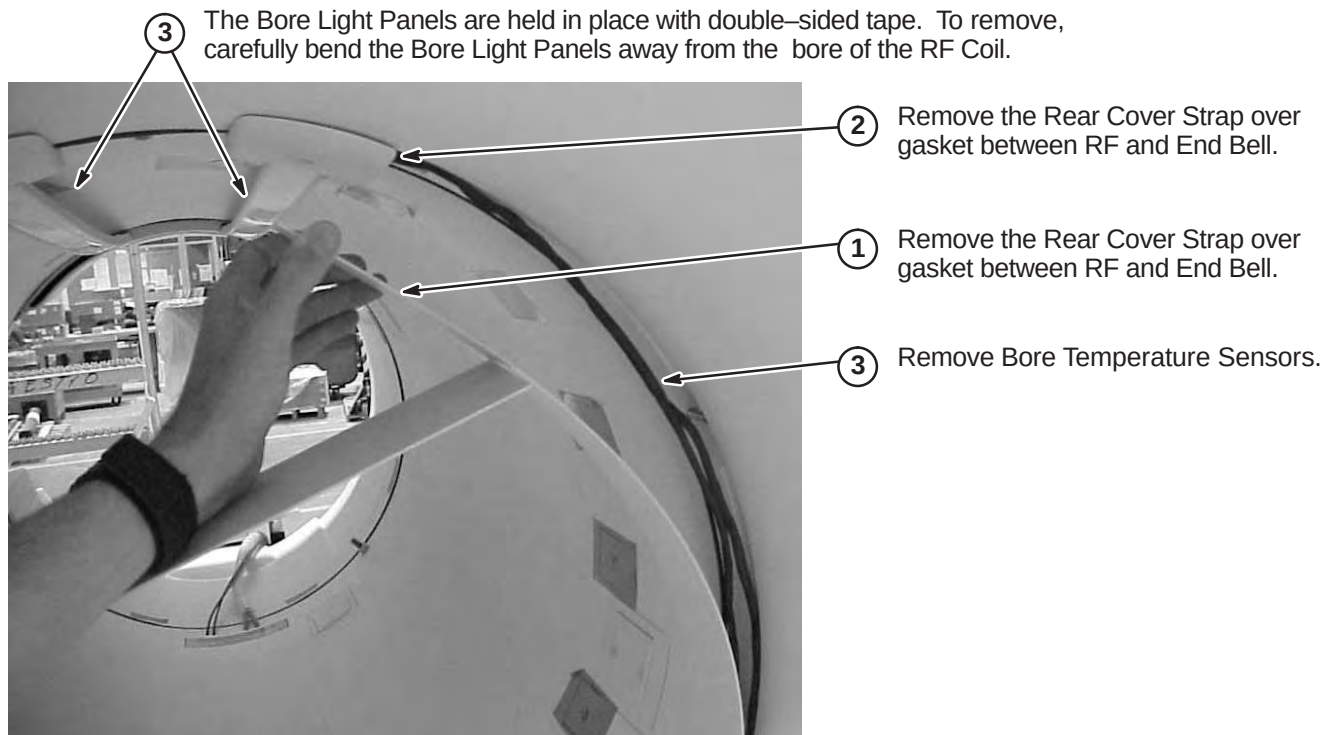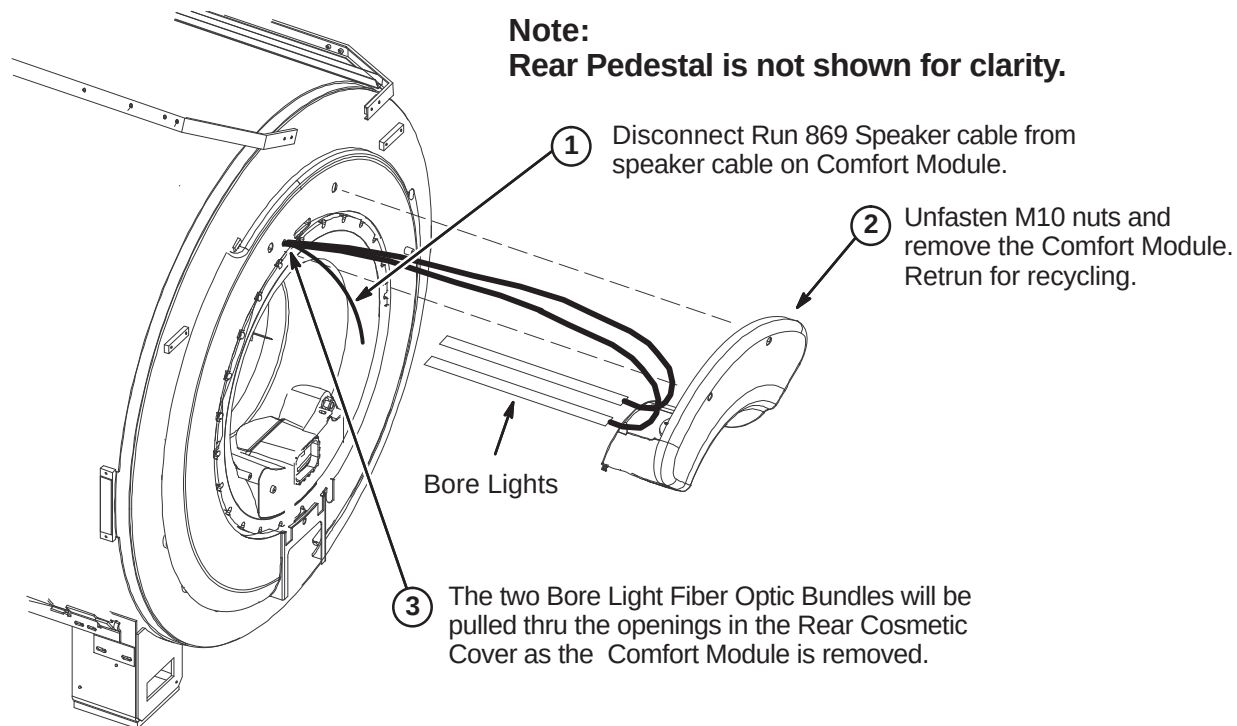4-1-19 Remove Rear Cover Strap and Bore Light Panels**4-1-20 Remove Comfort Module**


4-1-21 Cut the Gradient Cables

WARNING!**FERROUS MATERIAL HAZARD!**

Gradient cable cutters contain ferrous materials. Use extreme caution when handling the gradient cable cutters in the magnet room. Two people must be present at all times and tool must be kept outside the 200G field and as far from the magnet as possible. Refer to Section 1-1.

CAUTION


Do not cut gradient cables attached to the TRM Body Coil. Cut only the Gradient Output cables routed from Penetration Panel.

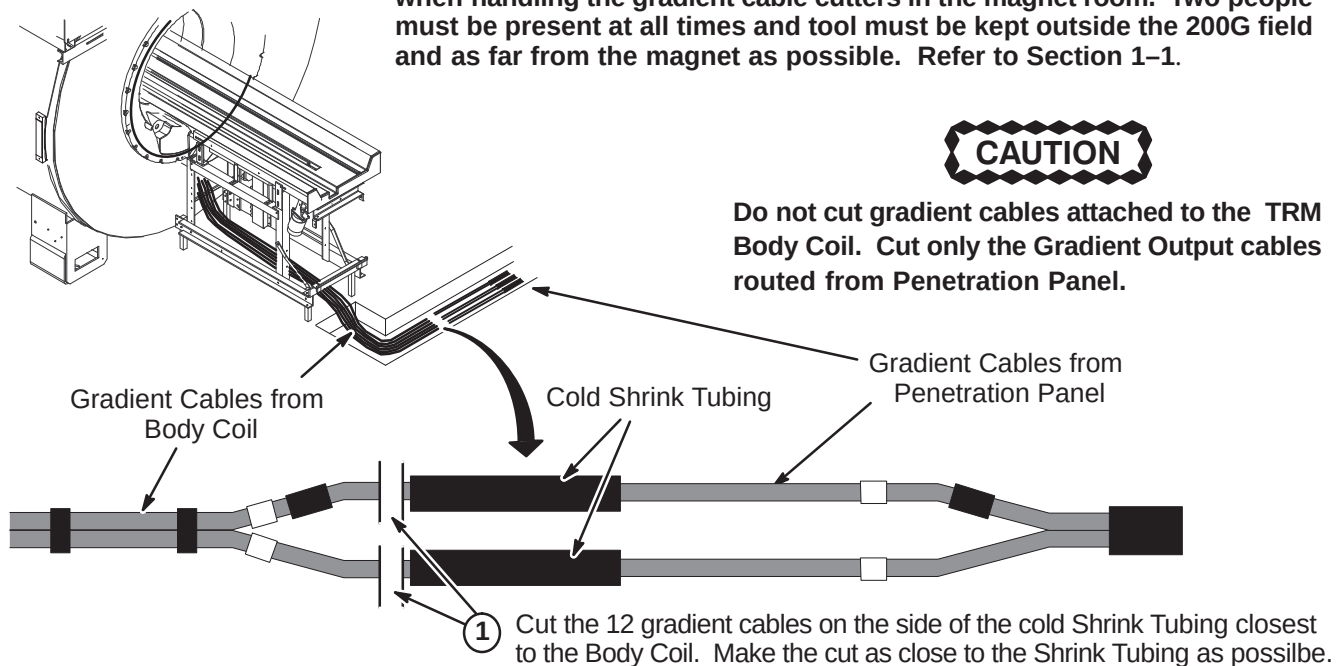

4-1-22 Remove Rear Pedestal Mounting Bracket

NOTICE

Only showing Mounting Bracket on one side. Repeat process for other side.

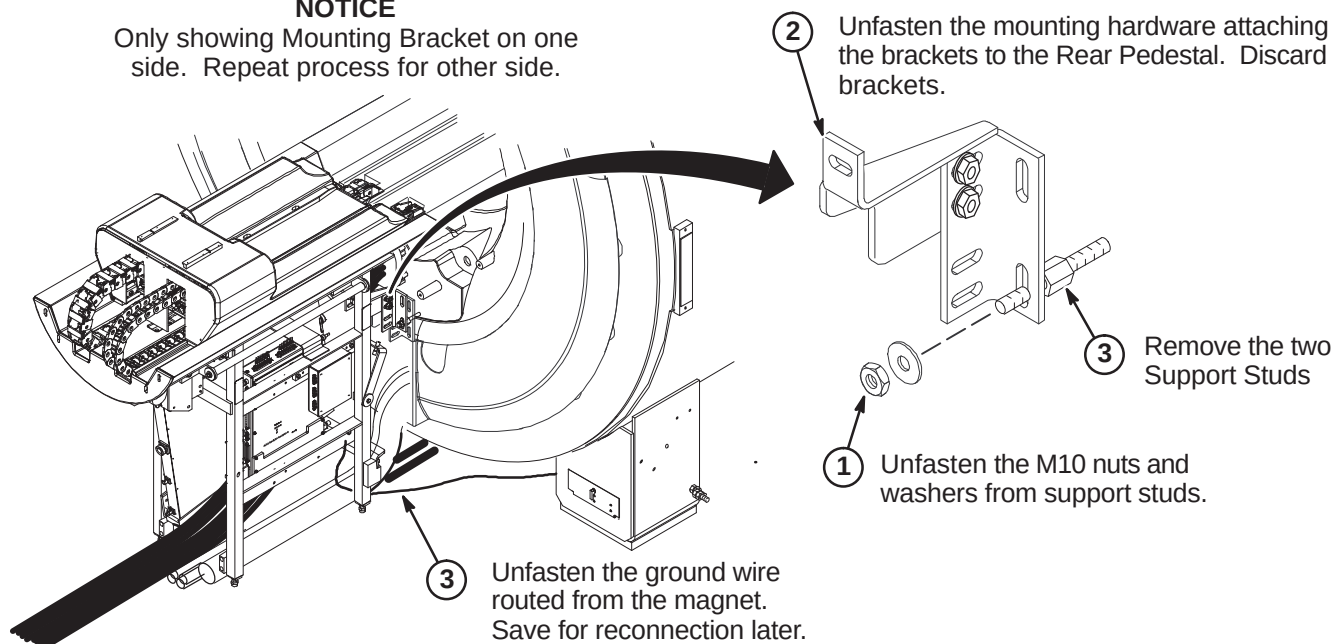

4-1-23 Move Rear Pedestal Away from Magnet

WARNING!**FERROUS MATERIAL HAZARD!**

REAR PEDESTAL REQUIRES AT LEAST THREE PEOPLE TO REMOVE FROM THE MAGNET ROOM. THE LONGITUDINAL DRIVE MOTOR IS FERROUS. KEEP THE REAR PEDESTAL, ESPECIALLY THE LONGITUDINAL DRIVE MOTOR, AS FAR FROM THE MAGNET AS POSSIBLE AT ALL TIMES.


CAUTION

Use care to not damage modules in Rear Pedestal when sliding the Rear Pedestal past the gradient cable clamp blocks.

- ① Move Rear Pedestal away from magnet to allow room for removing Rear End Bell.


Note:

For ease of moving Rear Pedestal, more cables routed from Penetration Panel may need to be disconnected. These cables will be reconnected during the upgrade process.

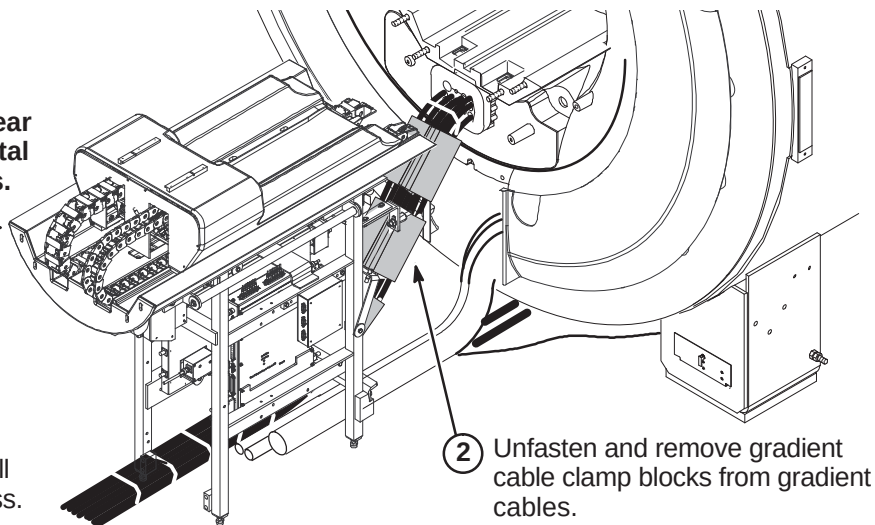


- ② Unfasten and remove gradient cable clamp blocks from gradient cables.

4-1-24 Remove Rear Cosmetic Cover

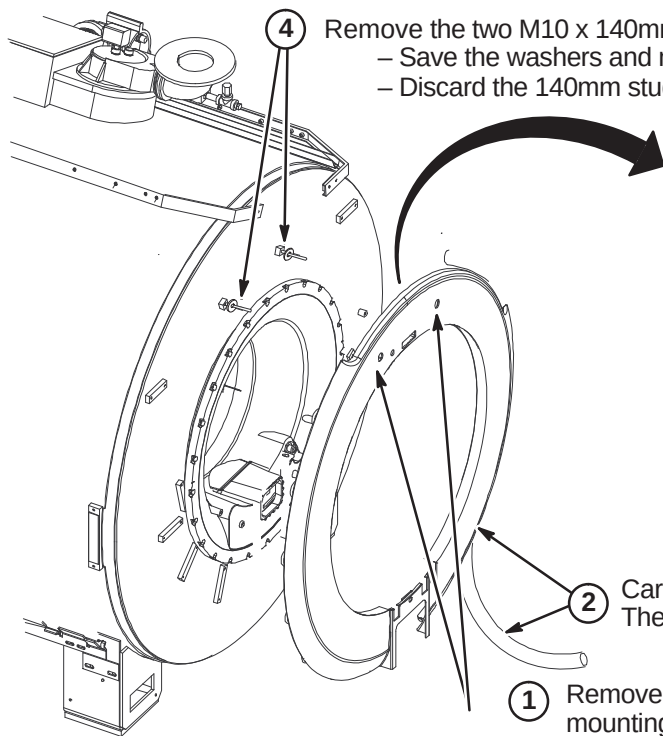

- ④ Remove the two M10 x 140mm Studs, washers, and nuts.
 - Save the washers and nuts for reuse
 - Discard the 140mm stud. A shorter stud will be installed for the upgrade.

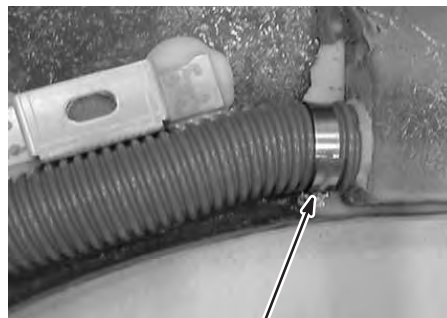





- ③ On back side of cover, unfasten the hose clamp and disconnect the 2 inch dia hose. Retain hose for re-installation later in upgrade.
- ② Carefully remove Rear Cosmetic Cover from studs. The 2 inch diameter hose is still attached.
- ① Remove Nylon Washer and M10 Locking Nut on each mounting stud. Save for re-use.

4-1-25 Remove Rear End Bell Pang Block Gaskets


① Loosen all the nuts and set screws that were required for centering the RF Coil within the TRM Coil. **DO NOT REMOVE** at this time.


② Unfasten the (18) M6 Pan Head Screws around the Block Cover

③ Unfasten the (16) M4 x 25mm Pan Head Screws in the inner holes of the Block Cover.

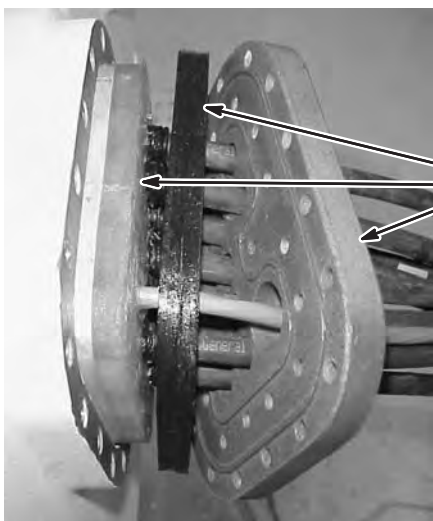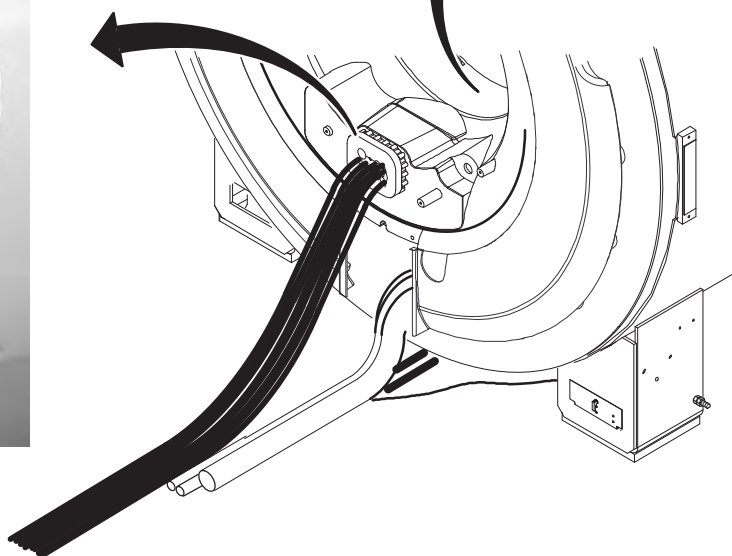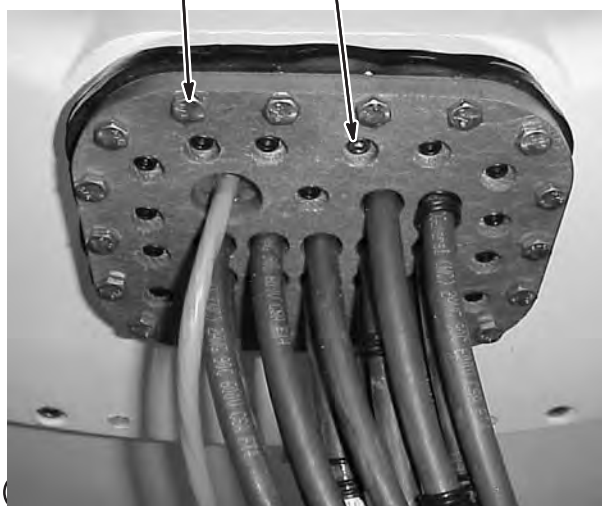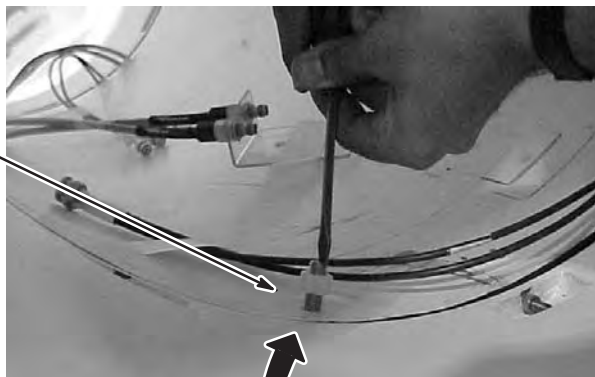



④ Remove all three gaskets by sliding them over the ends of the gradient cables. Return for recycling.

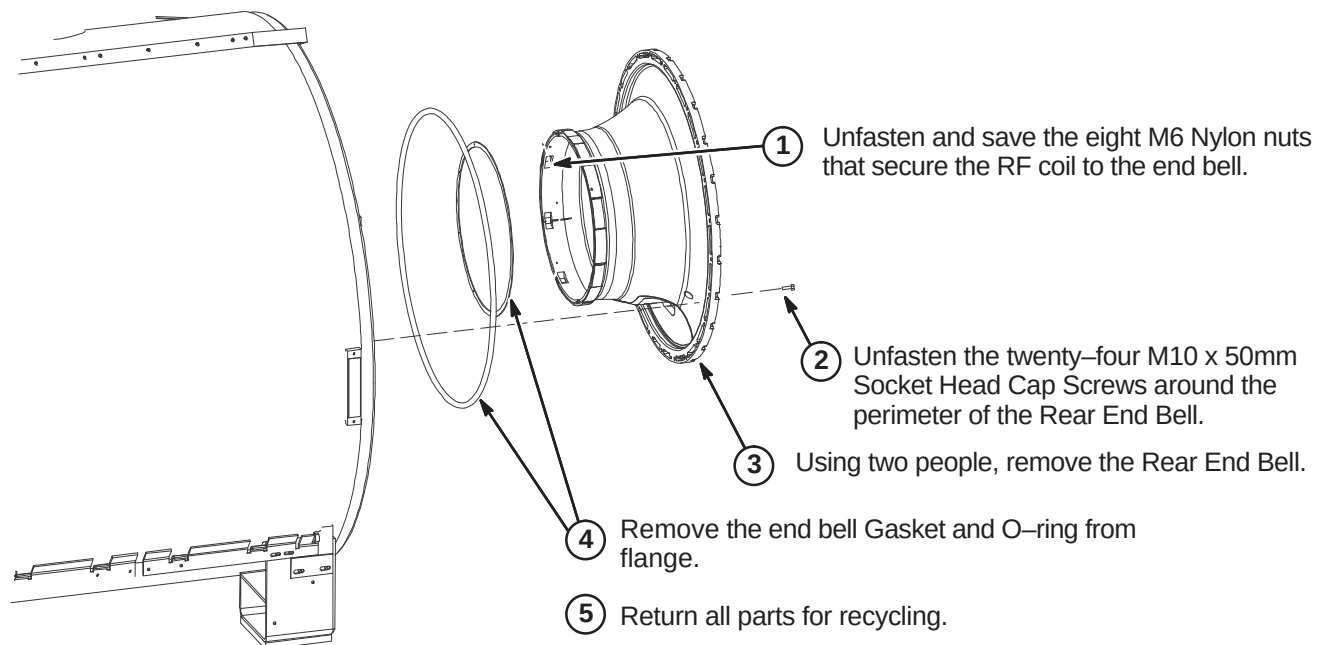4-1-26 Remove Rear End Bell and Gaskets


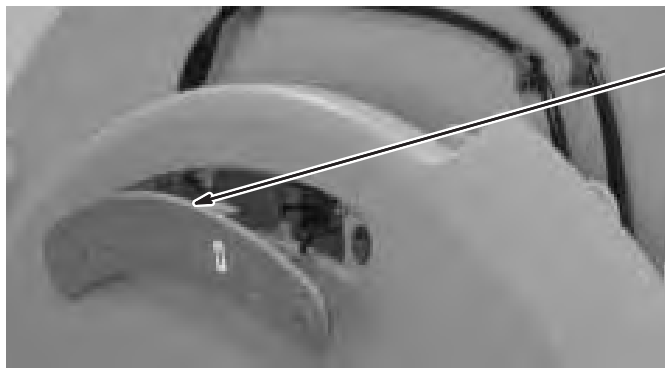
4-1-27 Remove Front Cosmetic Cover

- ① Carefully, using a small screwdriver, pry the Operator Display away from the front cover.
- ② Disconnect Run 1261 from J4 on back of Operator Display. Disconnect cables connecting to left and right Control Panels
- ③ Remove Operator Display. Return with Cosmetic cover for recycling.

- ④ Unfasten M10 Nuts and Nylon Washers from mounting studs. Save for reuse.
- ⑤ Carefully remove the Front Cosmetic Cover from the mounting studs.
- ⑥ Disconnect Run 385. The new Front End Bell will also come with a Run 385. Site FE must determine whether to keep existing cable or remove it for installation of new cable.
- ⑦ Remove the M10 x 65mm mounting studs, washers, and nuts from magnet remove the Front Cosmetic Cover from the mounting studs.

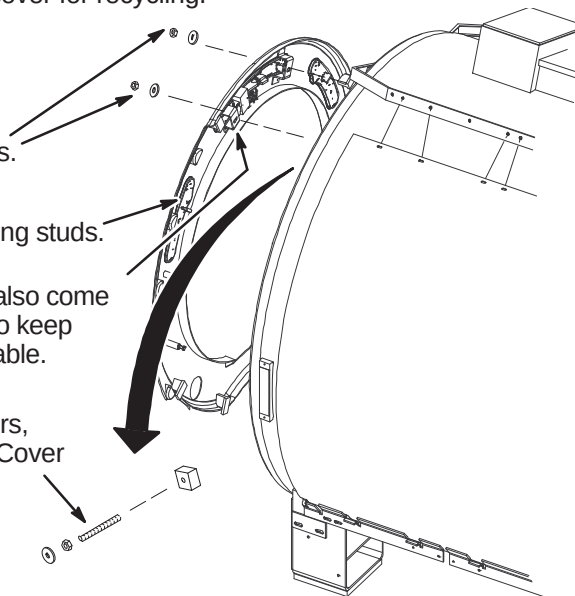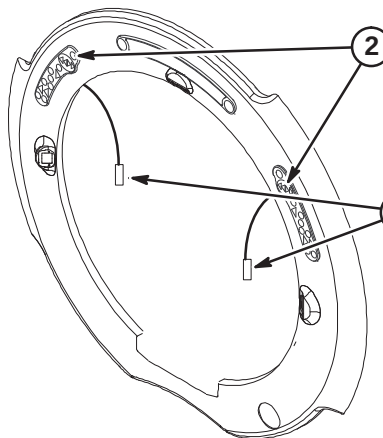**4-1-28 Remove Control Panels**

- ① Disconnect and save for transfer to new end bell the Left and Right Control Panel cables on Front Cosmetic Ring.
- ② Carefully, using a small screwdriver, pry the Left (2215716-3) and Right (2215716-2) Control Panels away from the front cover. Save for transfer to new front end bell.

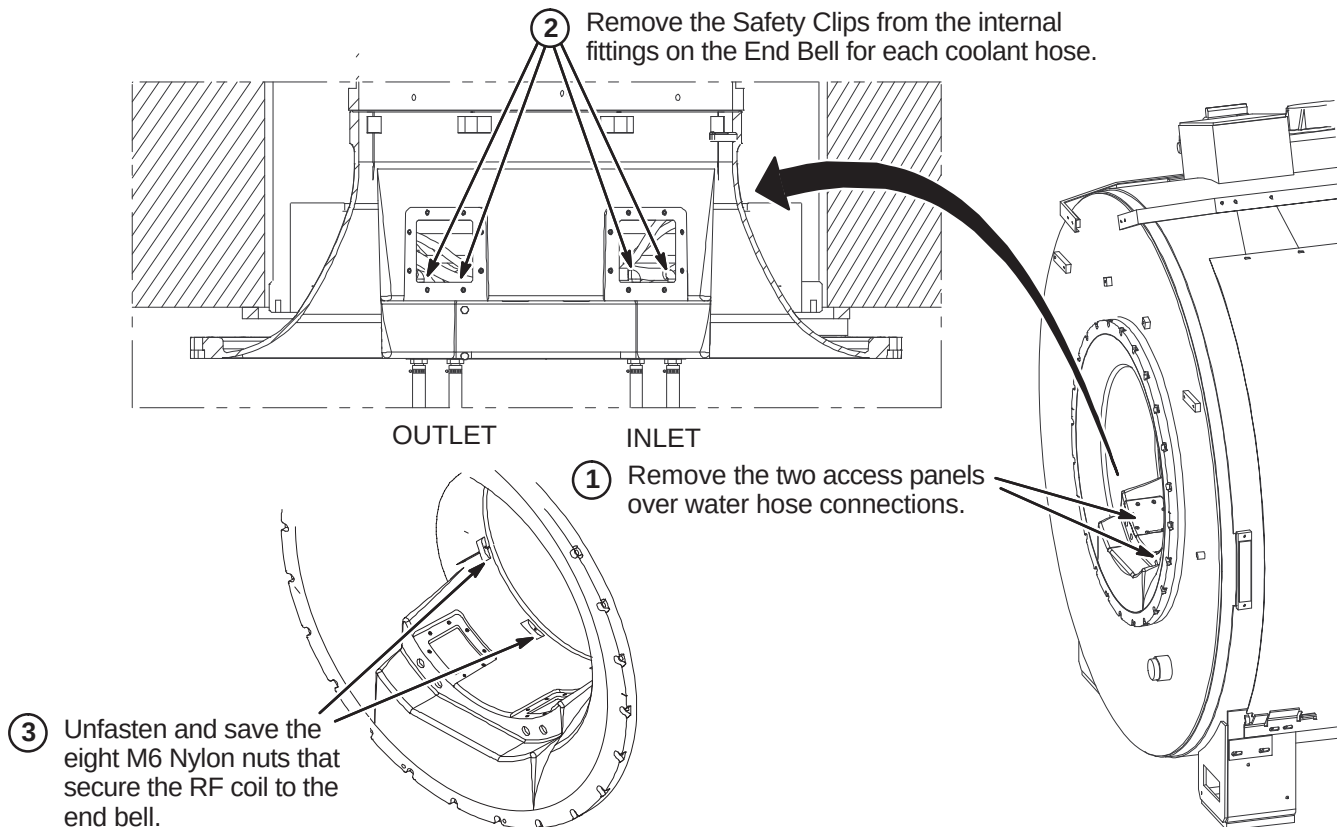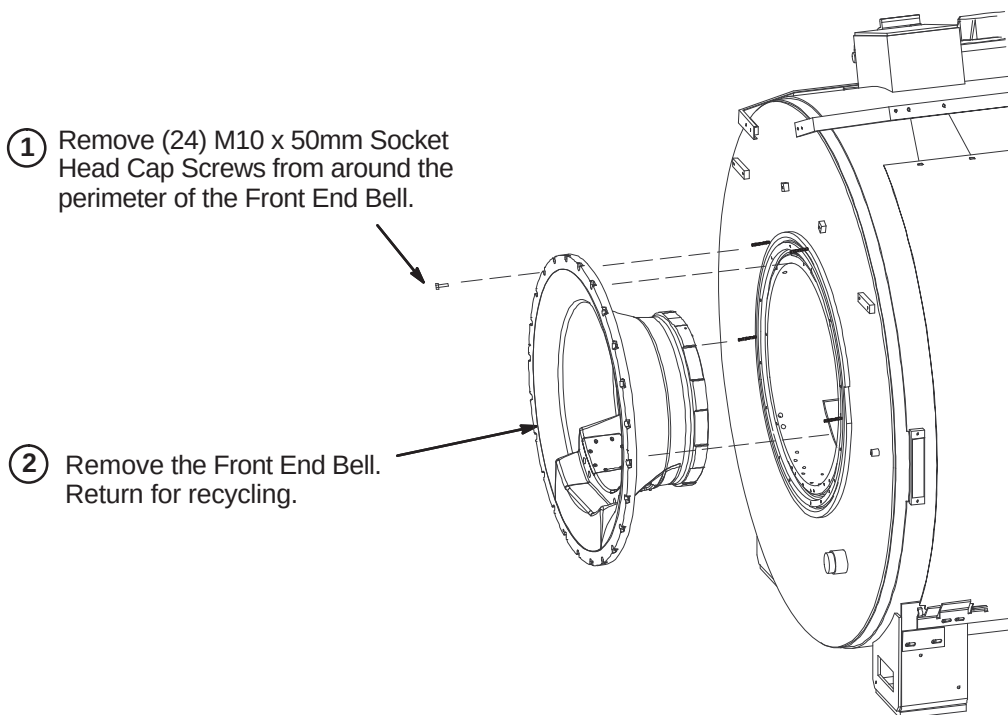4-1-29 Front End Bell Removal Preparation**4-1-30 Remove Front End Bell**


4-1-31 Removal of Gradient Coil


To remove the TRM Gradient Coil, refer to replacement procedures located on the MR Service Methods CD at:

For HD: Replacements>Coil Replacements>Twin Resonator Module (TRM) Replacements>Removing the TRM

4-1-32 Remove Magnet Interface Rings and TRM Supports

1. Remove TRM rubber pads and shim spacers
2. Unfasten M10 hex head Screw and remove TRM Mounting Support.
3. Unfasten the M10 x 25mm Socket Head Cap Screws around the perimeter of the Magnet Interface Ring.
4. Remove the Magnet Interface Ring.
5. Remove the Wedge Seal and Ring Gasket.
6. Remove the Wedge Ring.
7. Remove and clean any vacuum grease that may be on the magnet surface.
8. Repeat above steps for rear of magnet.

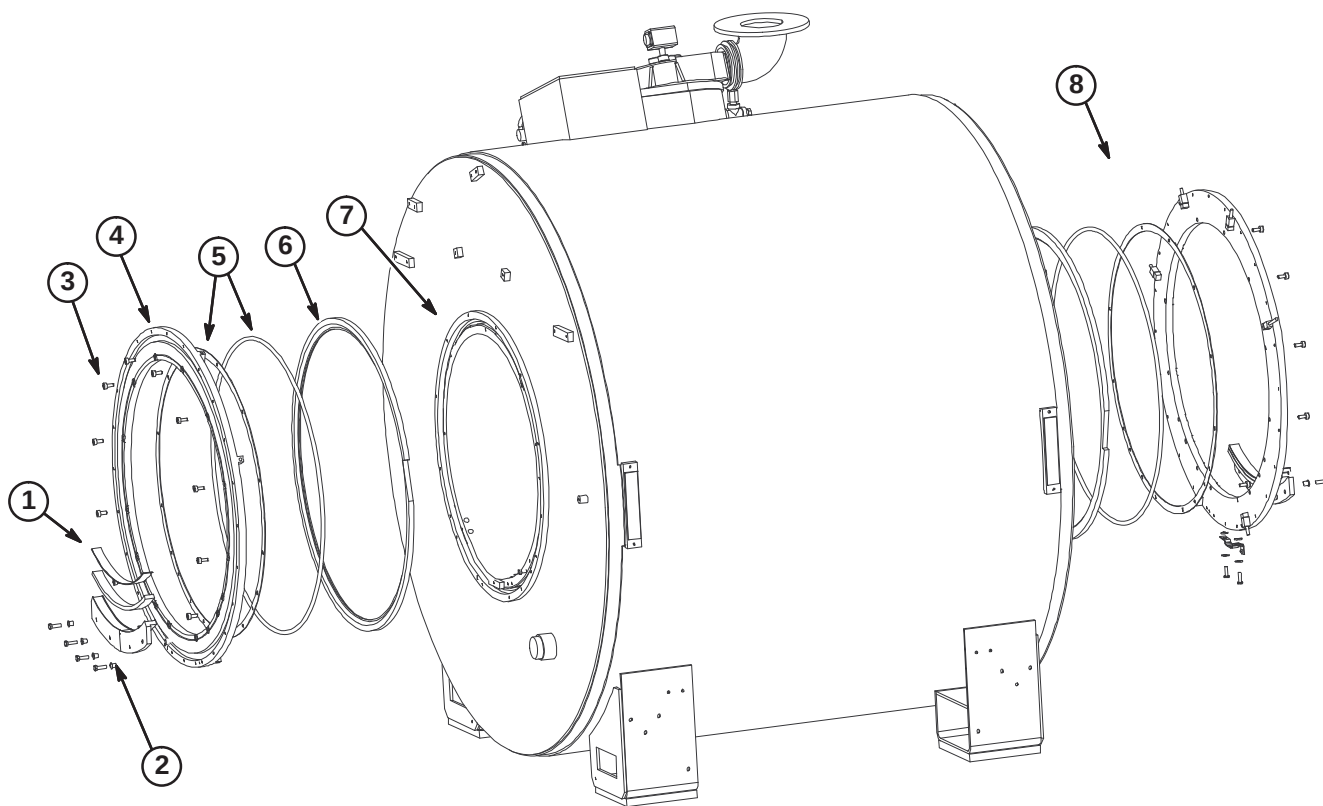

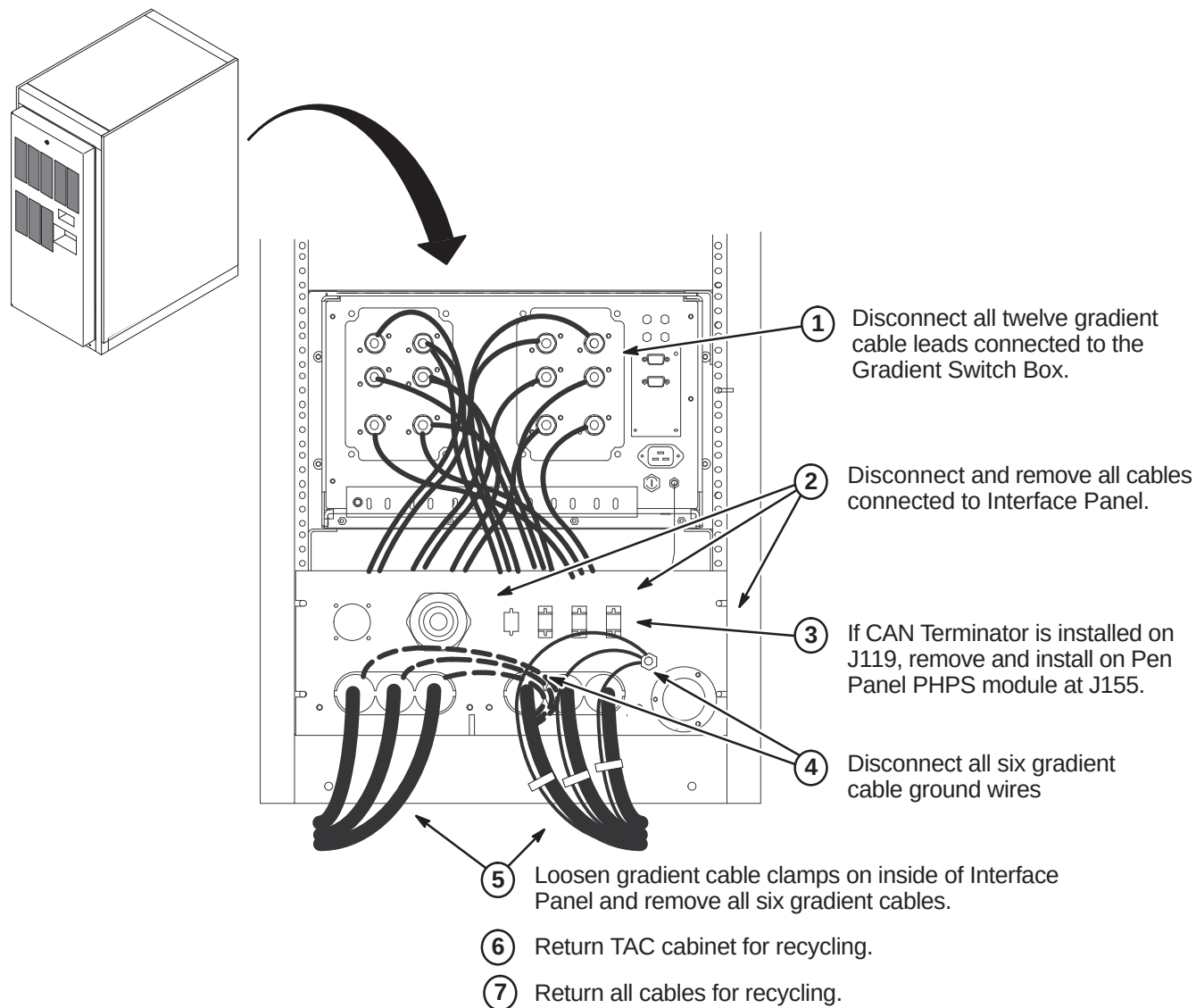4-2 DE-INSTALL TAC CABINET


4-3 DE-INSTALL GRADIENT CABLES



MAKE SURE ALL POWER TO CABINET IS OFF, LOCKED, & TAGGED BEFORE CONNECTING WIRES.


4-3-1 Disconnect Gradient Cables at HFD/PDU Cabinet

1. Unfasten and remove the top rear PDU module cover and ty-wrap to cabinet frame.
2. Unfasten the lower rear PDU interface panel, loosen cable clamps, and tilt top away from cabinet.
3. Loosen gradient cable clamps on inside of Interface Plate.
4. Loosen set screws securing gradient cable leads and remove gradient cables.

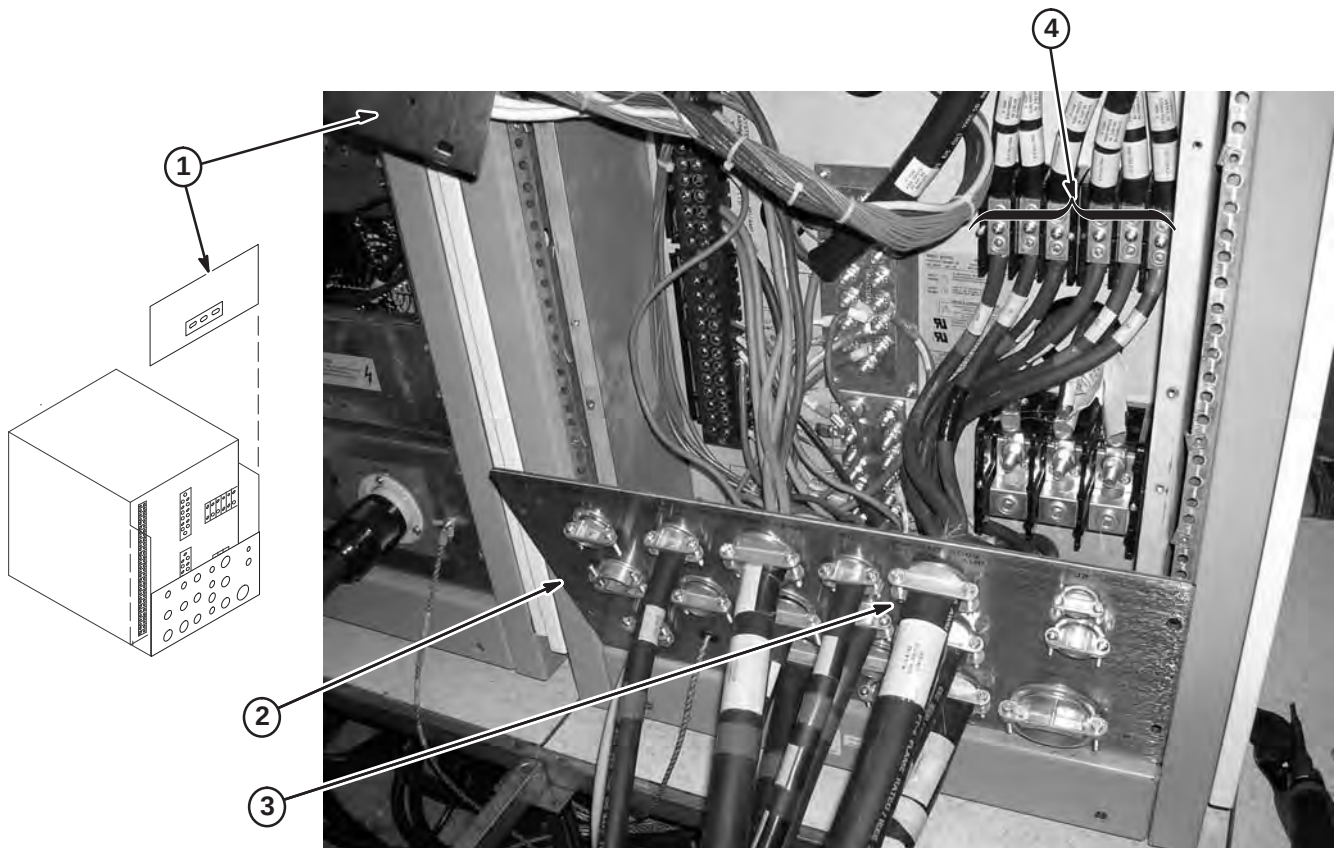

4-3-2 Disconnect Gradient Cables at Penetration Panel

1. Remove all gradient filter terminal covers on both sides of penetration panel.
2. Unfasten gradient cables from magnet room side of the panel.
3. Unfasten gradient cables from equipment room side of the panel.
4. Remove and return gradient cables for recycling.

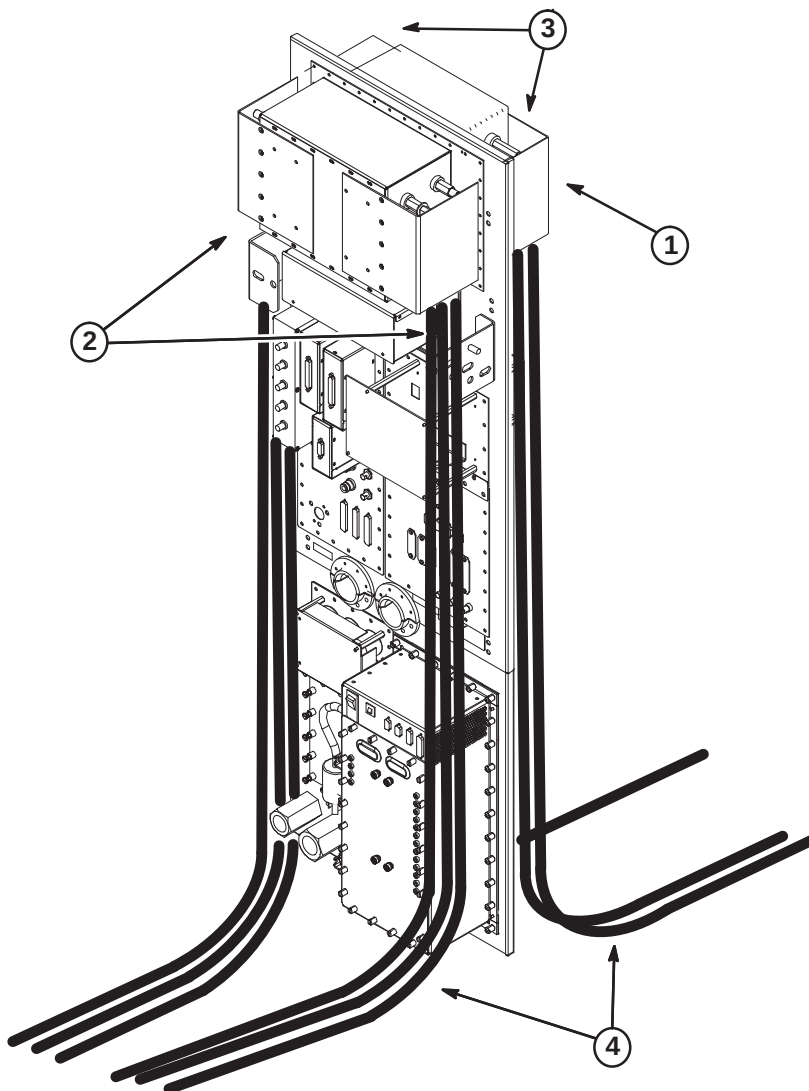


4-4 UPGRADE PENETRATION PANEL

1. Remove all gradient filter terminal covers on both sides of penetration panel.
2. Unfasten all screws attaching the TwinSpeed gradient filter to penetration panel. Save all fasteners. They will be used for attaching the BRM Gradient Filter.
3. **NOTE:** REFER TO WARNING BELOW. Remove gradient filter and return for recycling.
4. Install new Gradient Filter (46-320232G1) with 6-32 x 3/8 screws (46-208560P110) provided in the 5321590-4 upgrade kit.
5. If Shim Filter module was used for TwinSpeed, leave in place even though it will not be used with the upgrade.
6. If CAN terminator was installed at J119 on TAC, remove and connect to J155 on PHPS module.
7. Remove cover on Filter Panel on magnet room side. Connect Run 786 to F7 and F8 as marked. Route other end to Blower Box and connect to J2. Replace Filter Panel cover.

One of two styles on Penetration Panel may be installed on site. The older style without the PHPS module, or the newer style with the PHPS module. A different cable is installed for each style. Select either Step 8 or 9 to install.

8. **WITH PHPS:** Remove cover on Filter Panel on equipment room side. Connect Run 1183 to F7 and F8 as marked. Route other end to PDU in HFD cabinet and connect to Terminal Strip 30, Ground, and Neutral.
9. **WITHOUT PHPS:** Remove cover on Filter Panel on equipment room side. Connect Run 785 to F7 and F8 as marked. Route other end to System cabinet and connect to J47.
10. Replace Filter Panel cover on Penetration Panel.

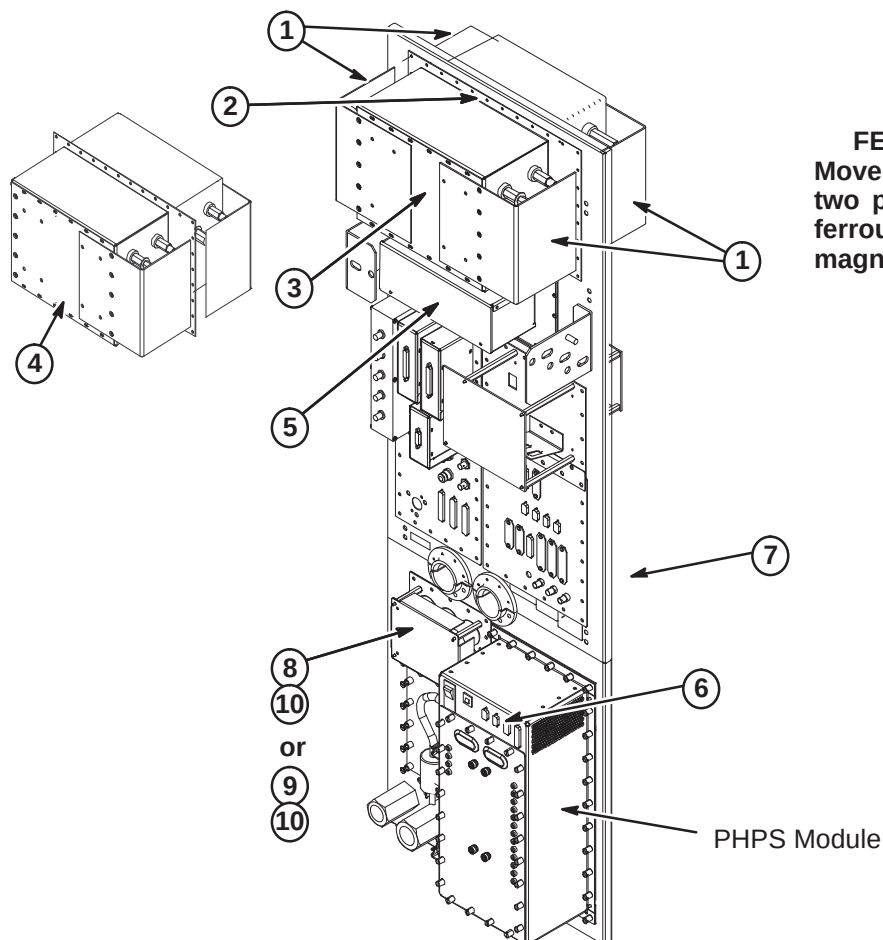
**WARNING!**

FERROUS MATERIAL HAZZARD!
Movement of the Gradient Filter requires two people. The Gradient Filter Box is ferrous and must be kept as far from the magnet as possible at all times.

4-5 INSTALL NEW MAGNET ENCLOSURE COMPONENTS


4-5-1 Install New BRM Body Coil

1. Install BRM mounting hardware to coil and magnet.
2. To install the new BRM Gradient Coil, refer to replacement procedures located on the HD MR Service Methods CD at:
Replacements>Coil Replacements>Gradient Coil Replacement for Cx Magnets>Gradient Coil Replacement Preparation and Gradient Coil Installation Procedure
3. Install Radial Supports between coil and magnet.

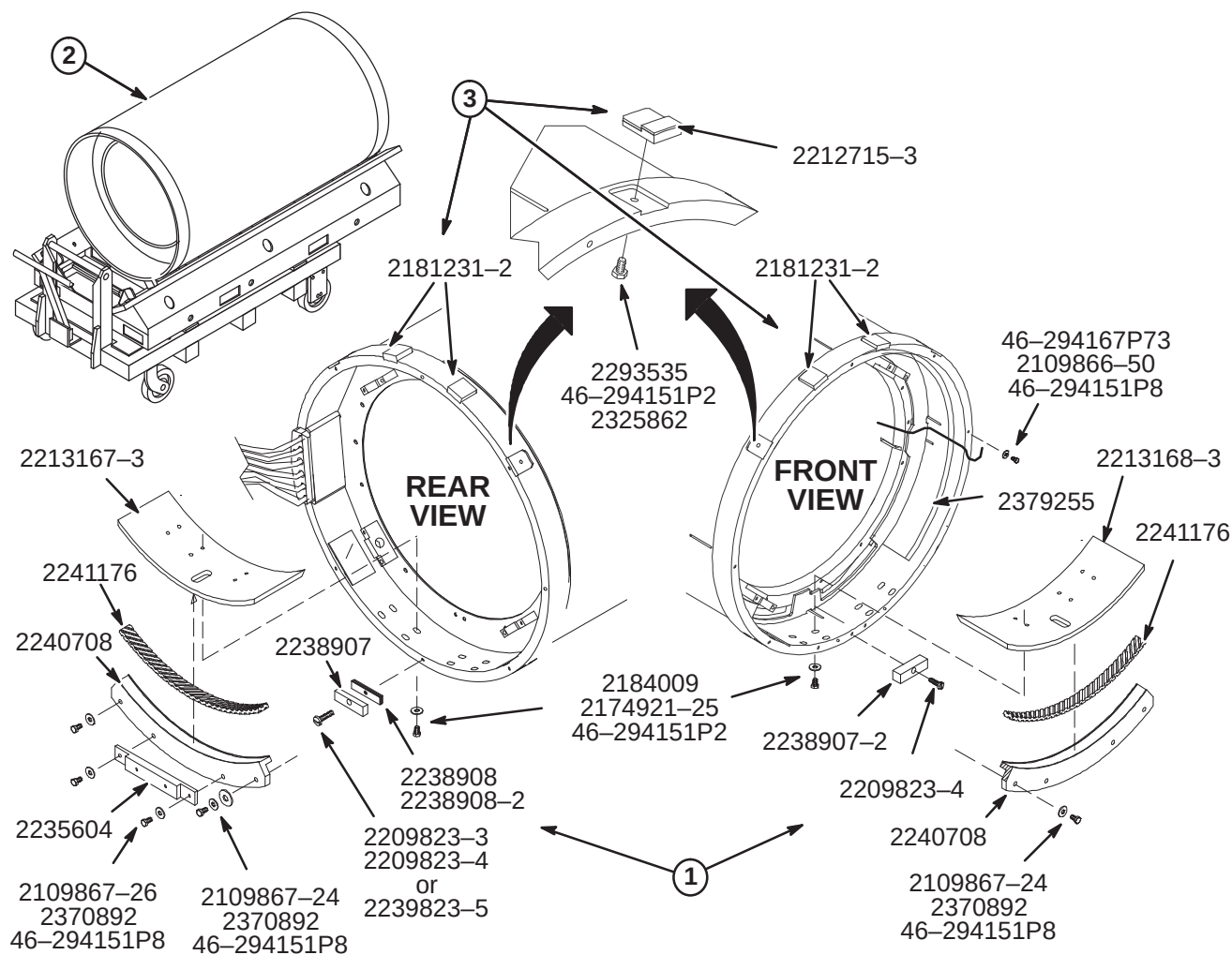


4-5-2 Install Gradient Cable Lead Board


Refer to Gradient Lead Board Kit, 5110336, shipped with catalog M3335BD, for list of parts and the install manual, Direction 2404642, *Signa Gradient Cable Lead Board Upgrade Kit*, for procedures to install the parts listed below.

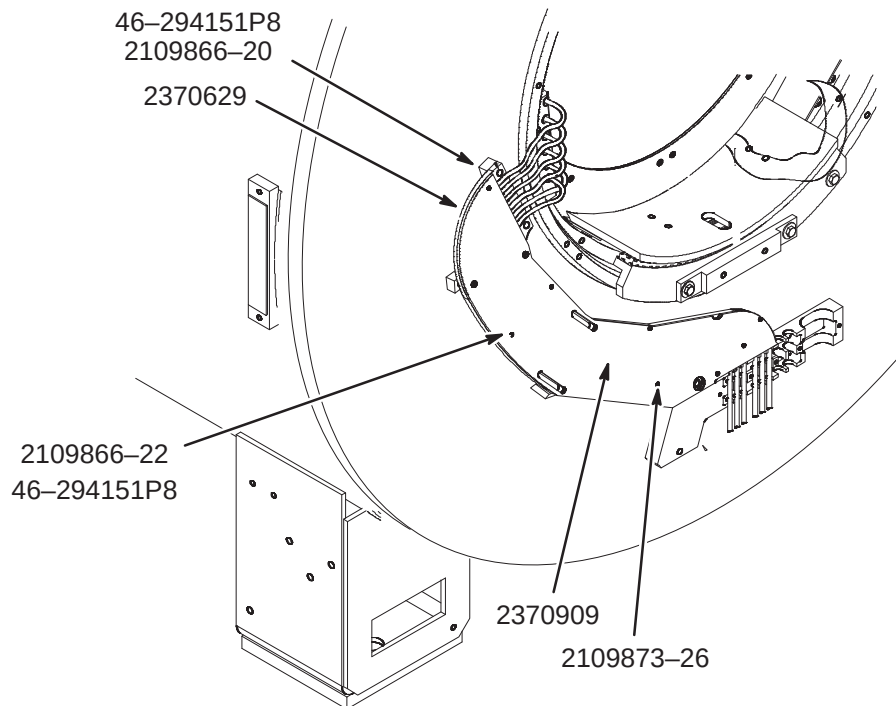


4-5-3 Install Cable Clamp Base


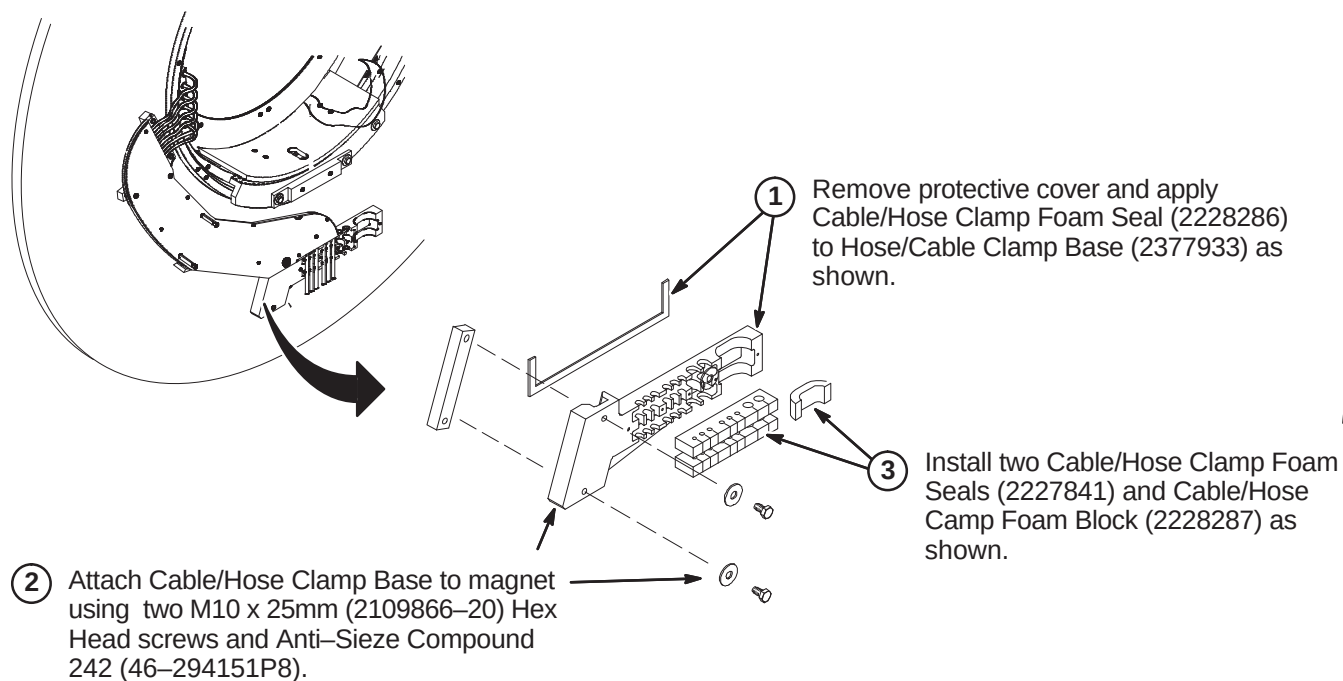




4-5-4 Connect Coolant Hoses to Gradient Coil

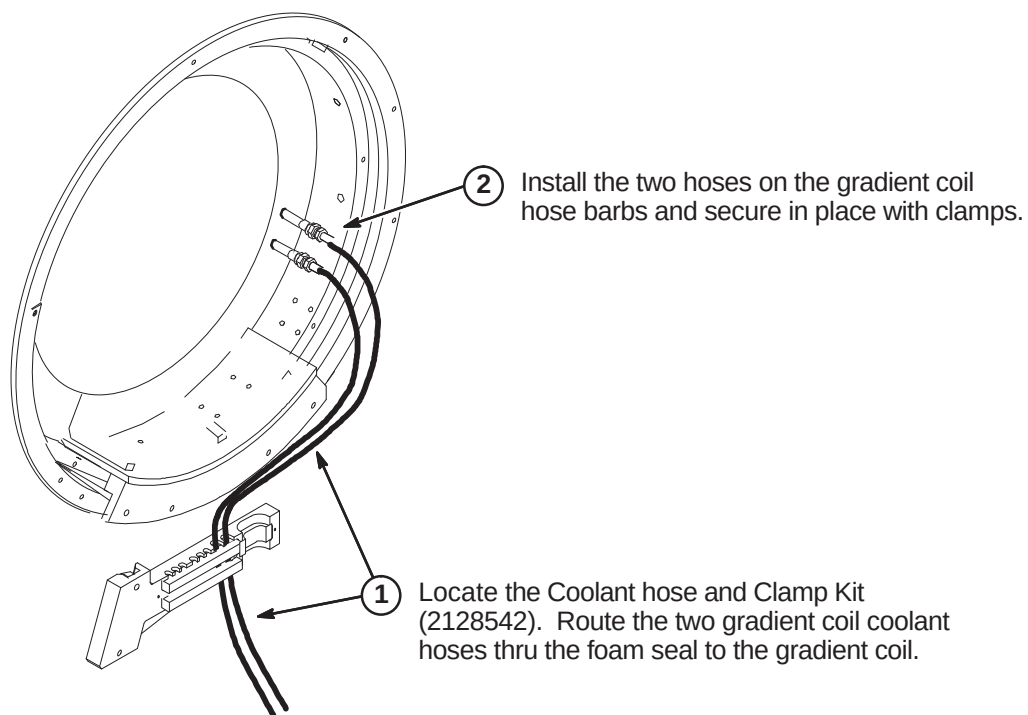

4-5-5 Install Cable/Hose Clamp

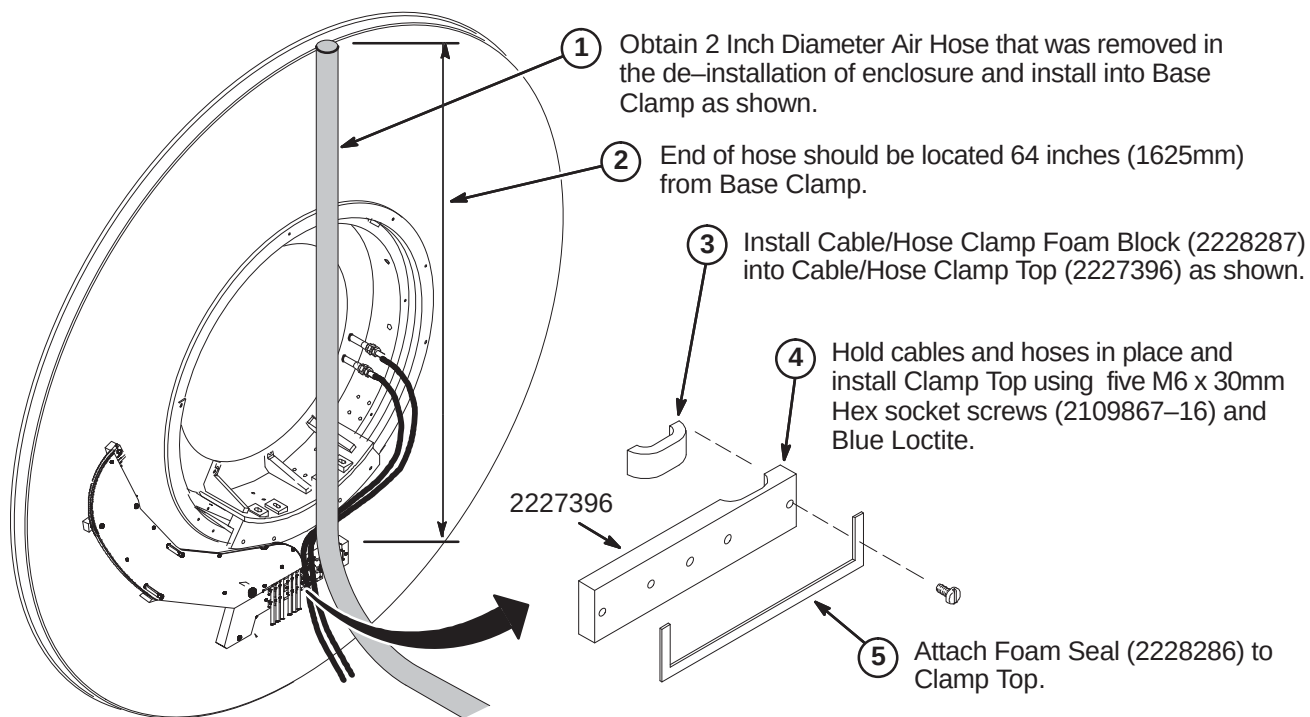


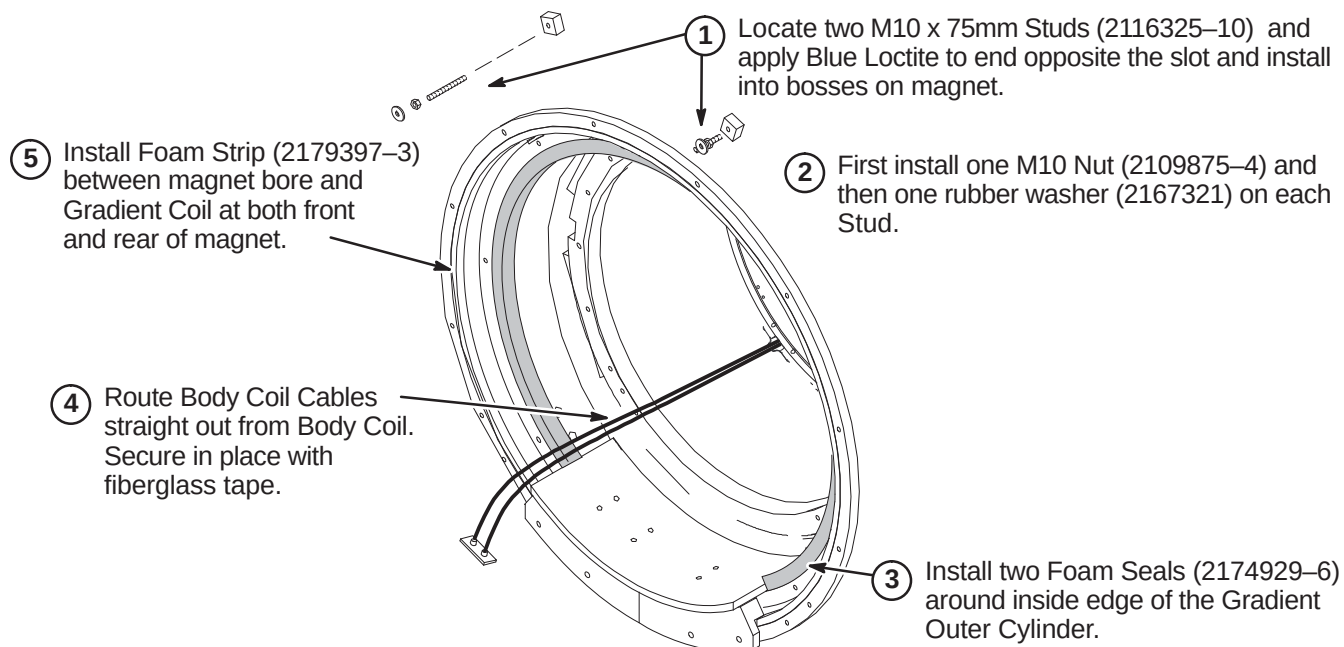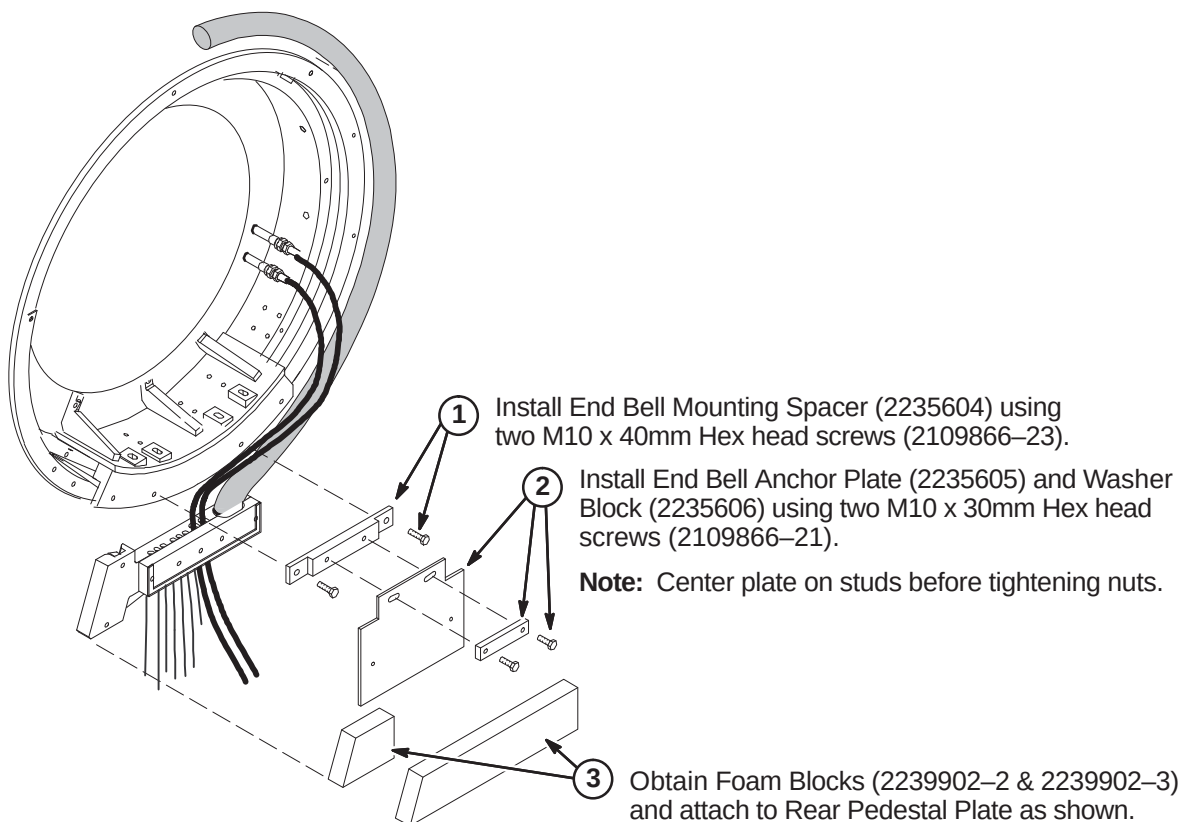4-5-6 Install Bore Front Seals, Mounting Bracket, and Studs**4-5-7 Install Foam Gasket and Gradient cable Air Seal**

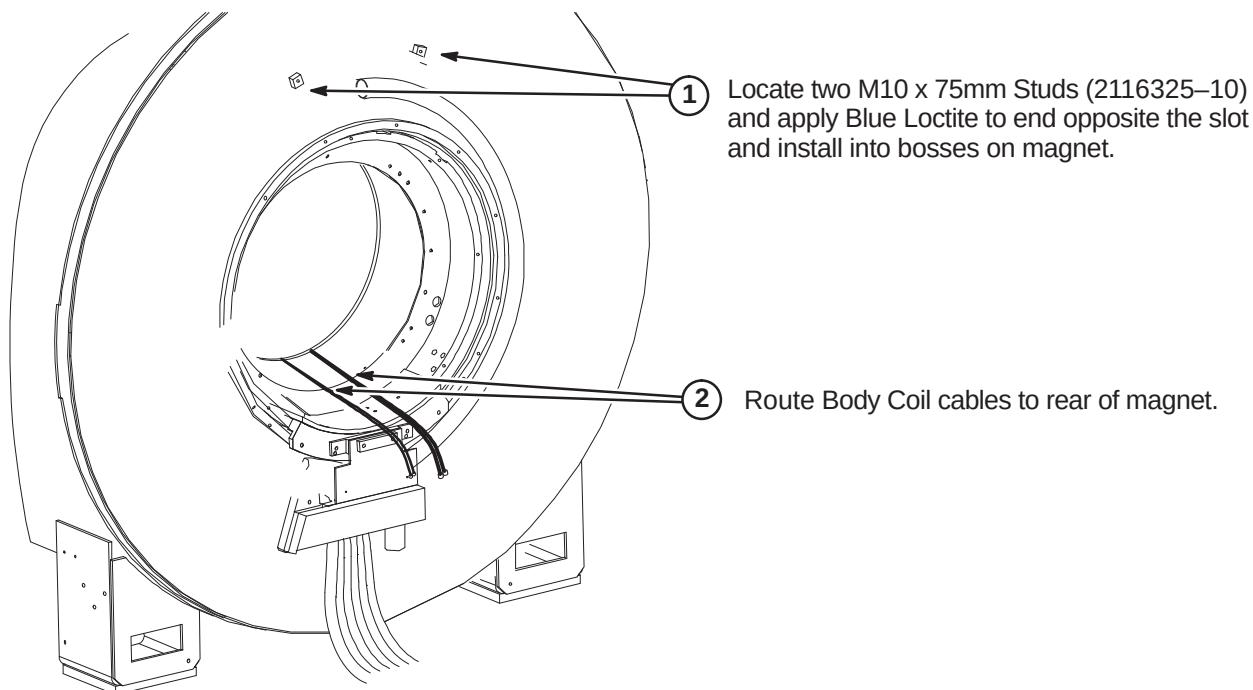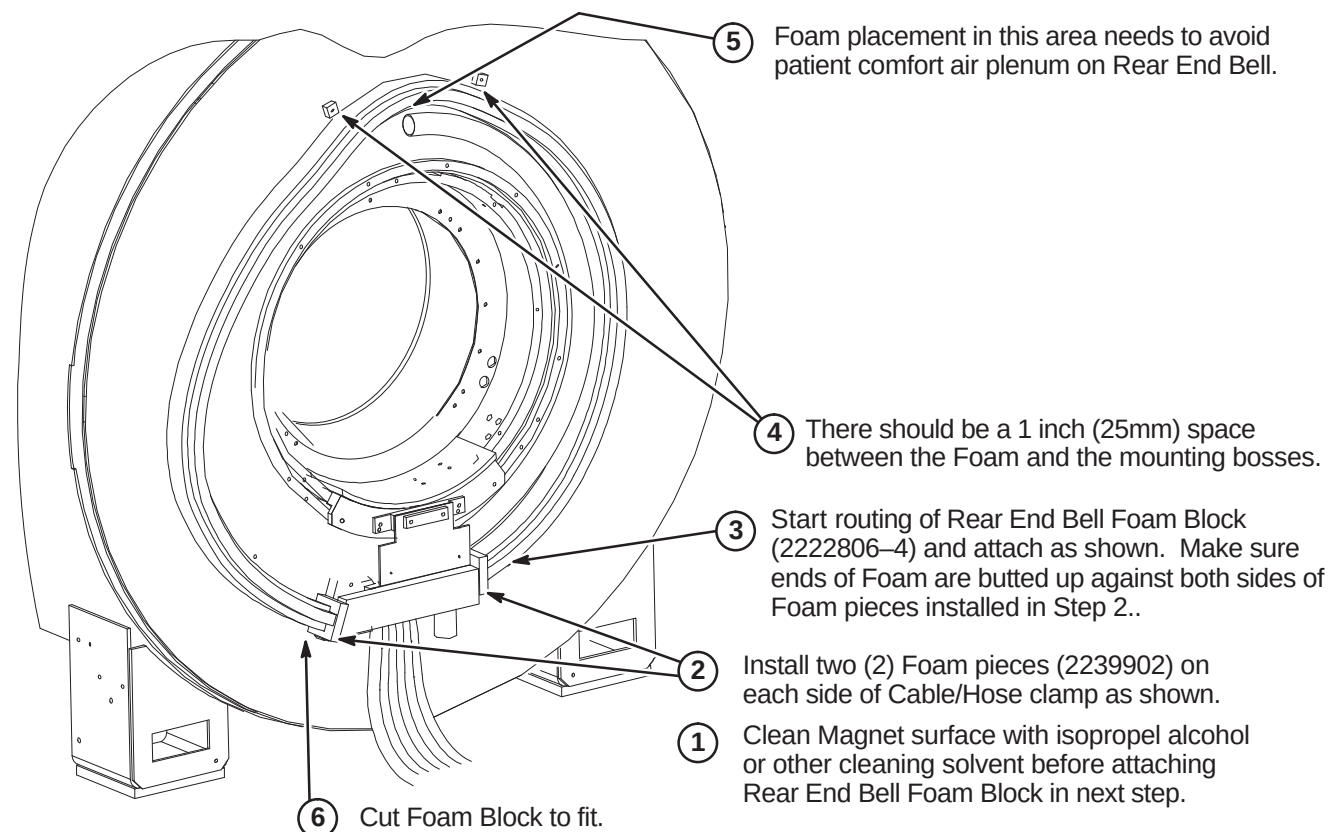4-5-8 Install Foam Gasket and Gradient cable Air Seal**4-5-9 Attach Foam Block**

4-5-10 Position Rear End Bell Near Rear of Magnet


- ① Position Rear End Bell (5120954) near rear of Magnet.
- ② Attach previously removed 2 inch Air Hose to the port at the top of the End Bell using clamp (46-208765P28).

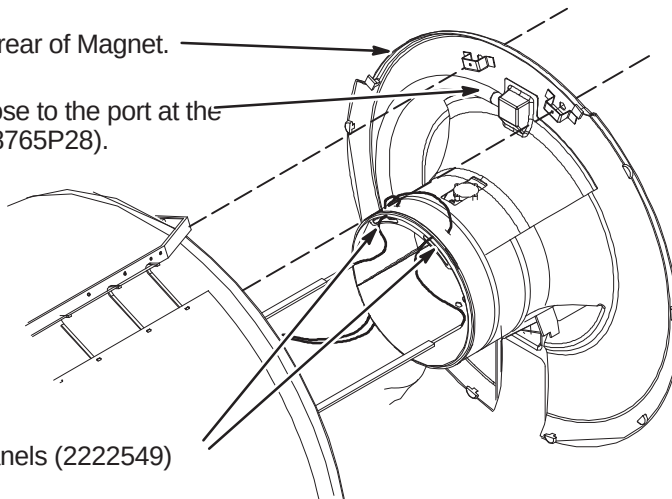


- ③ See Section 4-5-11, if Bore Light Panels (2222549) are not routed in End Bell.

4-5-11 Route Bore Light panel in End Bell


- ① Cut a small slit in two places of the thermal blanket for insertion of fiber optic cable (part of Bore Light Panel). Make sure adhesive side of Bore Light Panel is towards the top of RF Coil.
- ② Final installation will have Bore Light Panel inserted into position in Rear End Bell.
- ③ If possible, apply tape to slits to prevent escape of air flow.

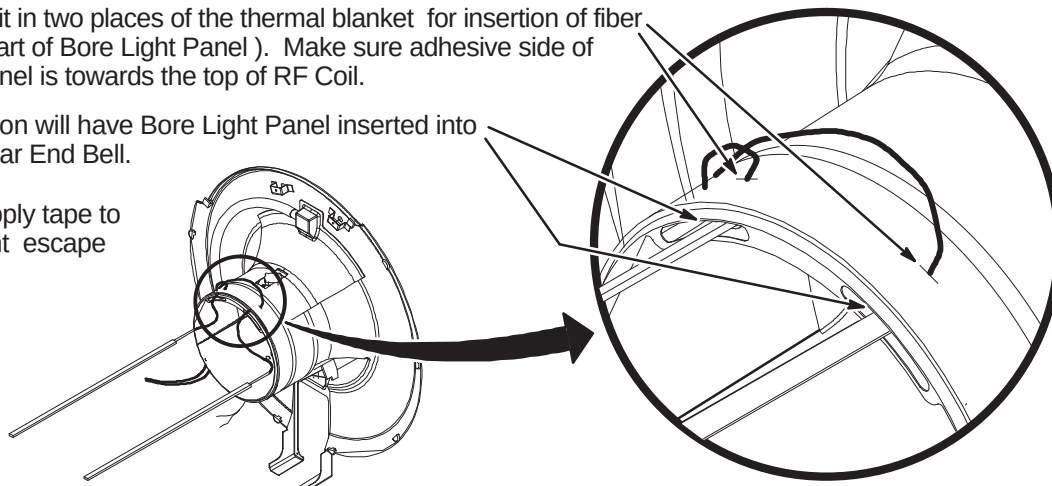**4-5-12 Route Body Cables thru End Bell Opening**

- ① Insert Body Coil cables thru rectangular opening in the bottom of EndBell.
- ② Continue to move End Bell into position while aligning hose, cables, and fiber optics.

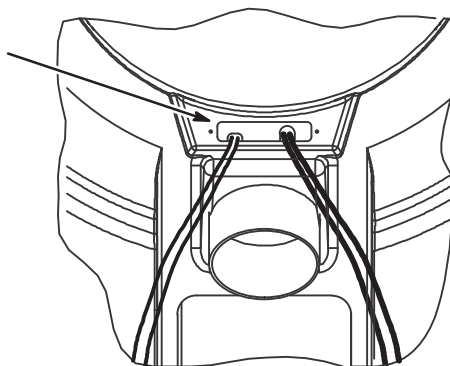



4-5-13 Secure Rear End Bell to RF Coil and Magnet

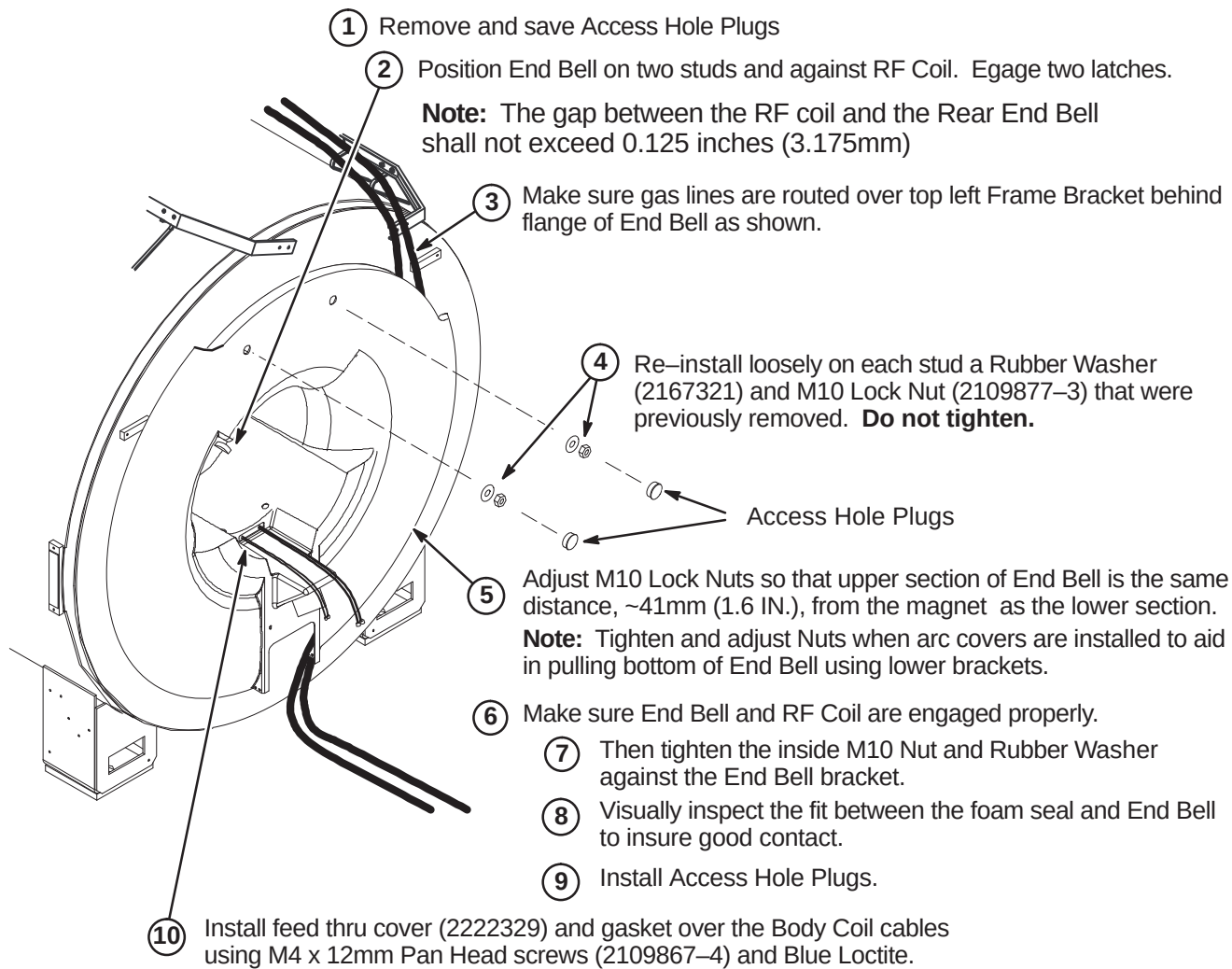

4-5-14 Install Bore Light Panels


- ① Remove the adhesive protector from the Bore Light Panels.
- ② Insert Bore Light Panels into Bore Light End Caps and tighten screws.
- ③ Press Bore Light Panels into place.

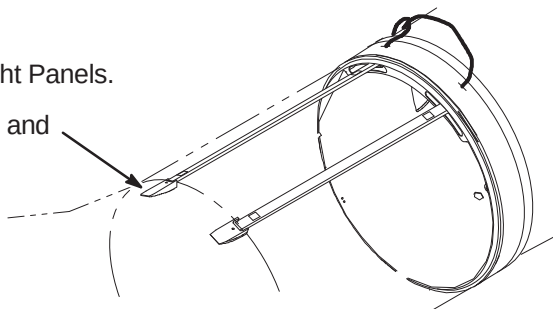

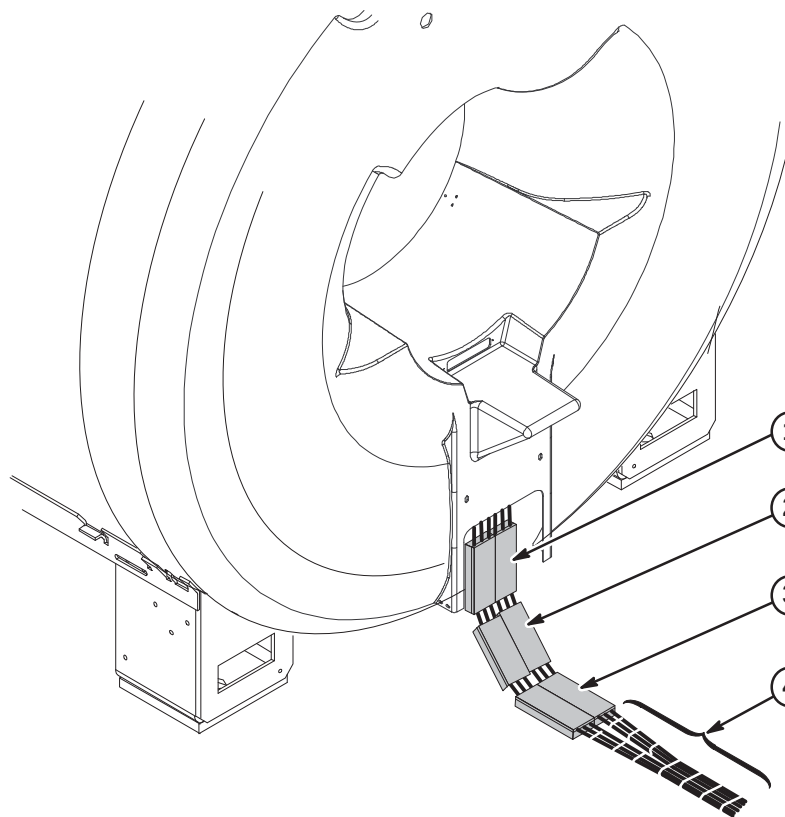4-5-15 Install Gradient Cable Clamp Blocks


Note: Locate 5182681-2, Gradient Cable Clamp Block kit. Six blocks (5179382) and six covers (5179303) will be used.

Orientate blocks as shown for installation steps below.

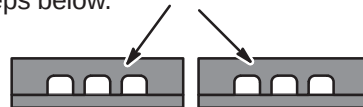

- ① Install top row 50mm to 75mm (2 to 3 inches) from the end of trough.
- ② Install middle row 25mm to 50mm (1 to 2 inches) from top row.
- ③ Install bottom row 25mm to 50mm (1 to 2 inches) from middle row.
- ④ Install tie-wraps, about 3 to 4 inches apart (75mm to 100mm) around the remaining lengths of the gradient cables.

4-5-16 Install Split Bridge Front Support

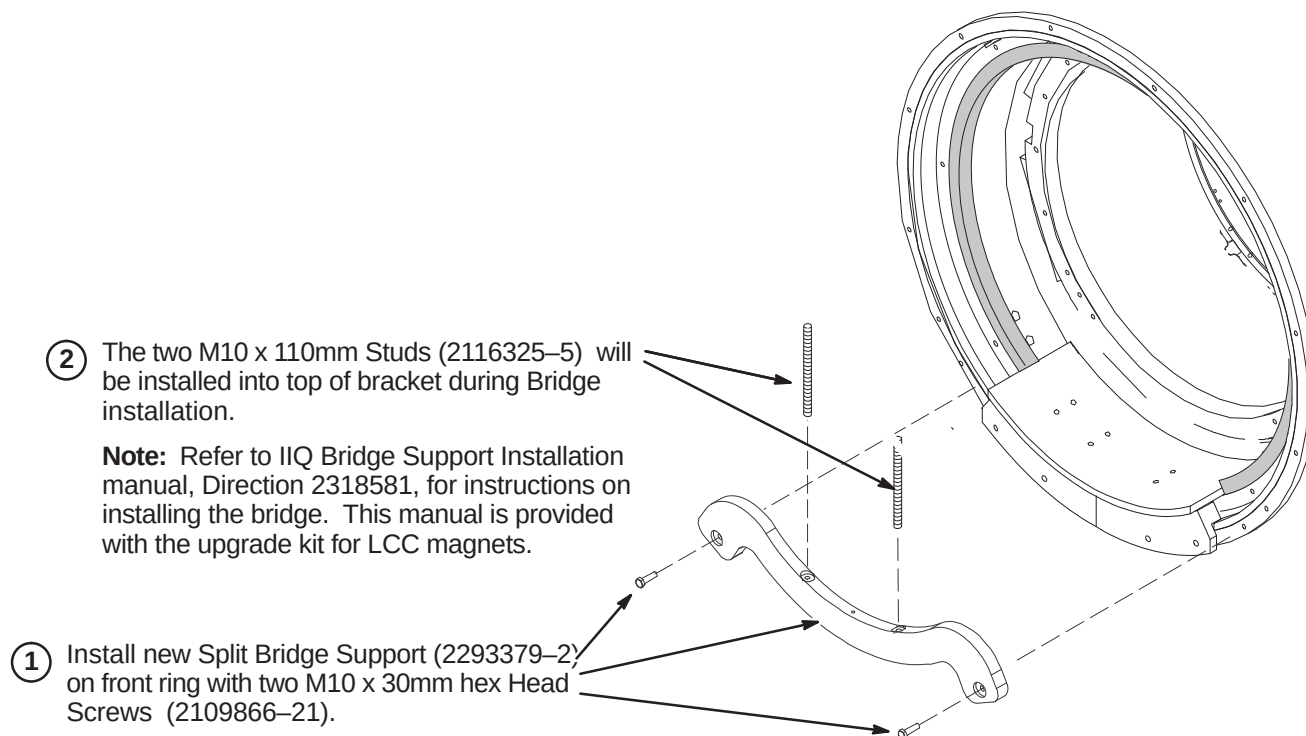


4-5-17 Position Front End Bell in Front of Magnet



- ① Carefully move the Front End Bell (2315088) into position in front of the magnet.
- ② Connect existing Runs 778 and 779 to RF Coil cables.
- ③ Attach Bore Light End Caps (2218899) to End Bell using two M4 x 12mm Hex Head Screws (2109867-4) and Blue Loctite in each End Cap.
- ④ Connect existing Run 385 to JP1 on Voltage Regulator board located behind Display Panel.



NOTE:


Make sure new Front End Bell has Dual Temperature Sensor cable, Run 1260, attached before installation.

- ⑤ Carefully pull Operator Display Panel away from mounting studs. Display panel is connected to cables from left and right Operator Control panels. Let Operator Display panel hang below opening.

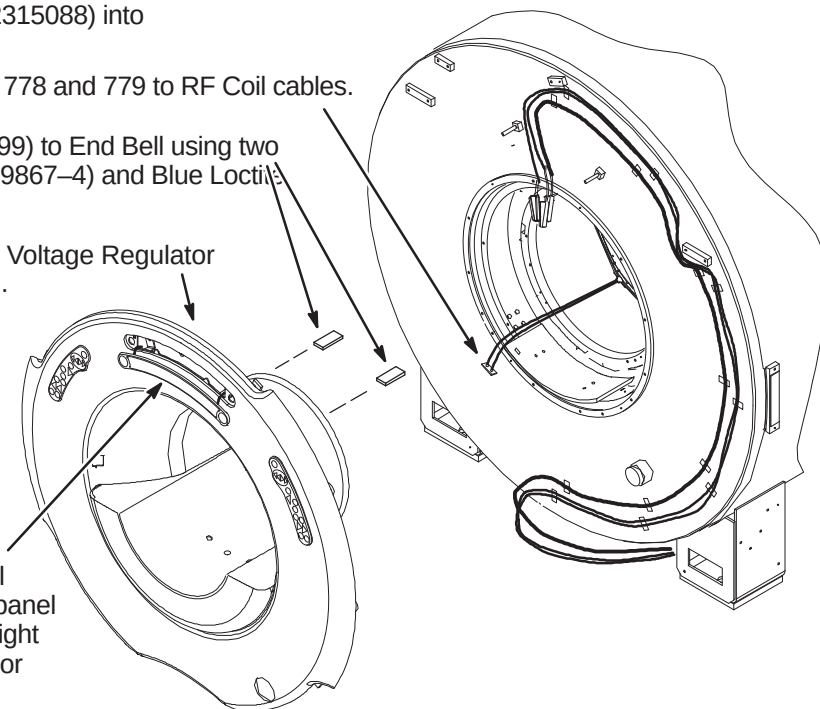

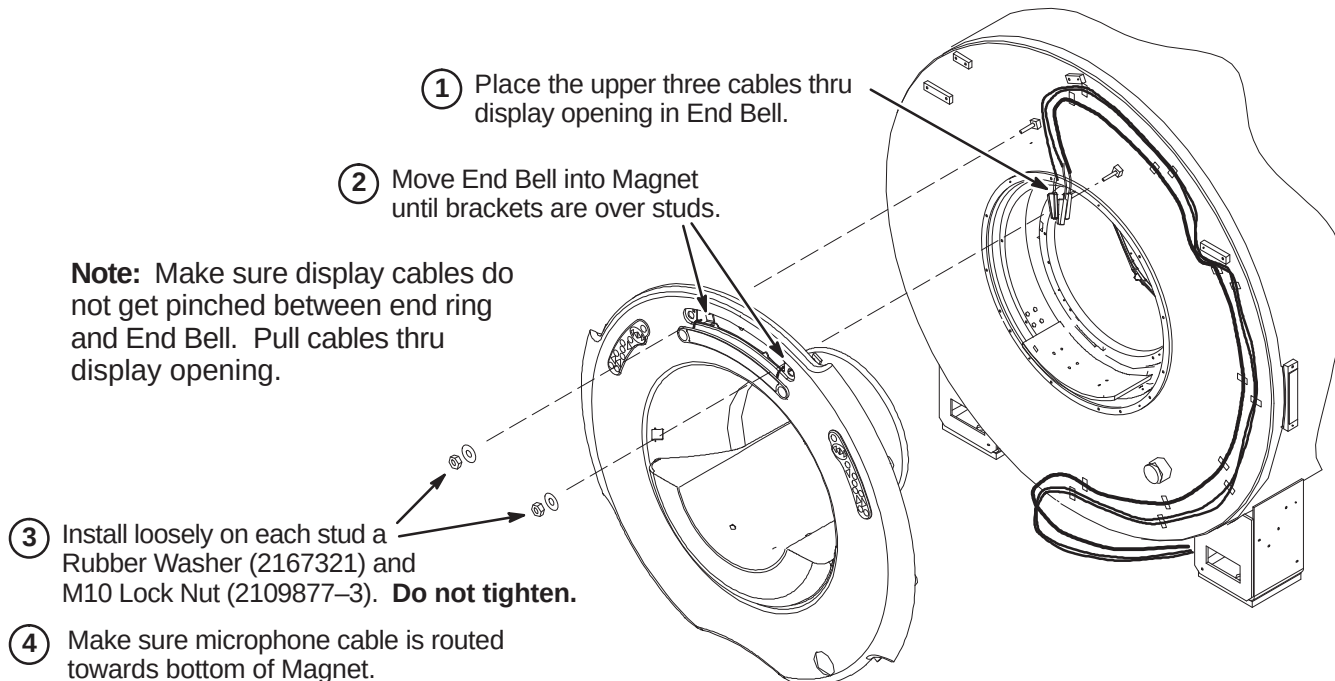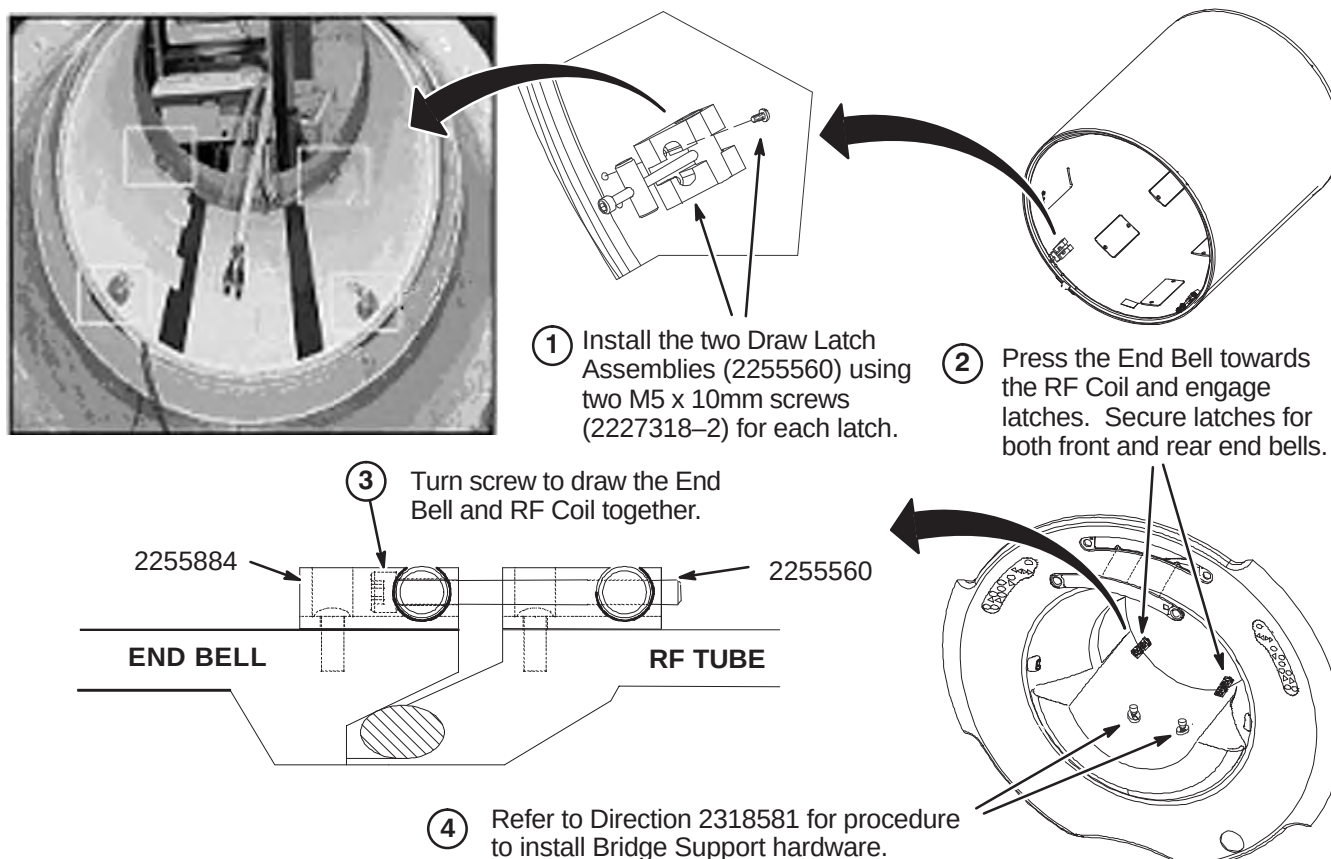4-5-18 Install Front End Bell on Magnet**4-5-19 Secure Front End Bell to RF Coil**

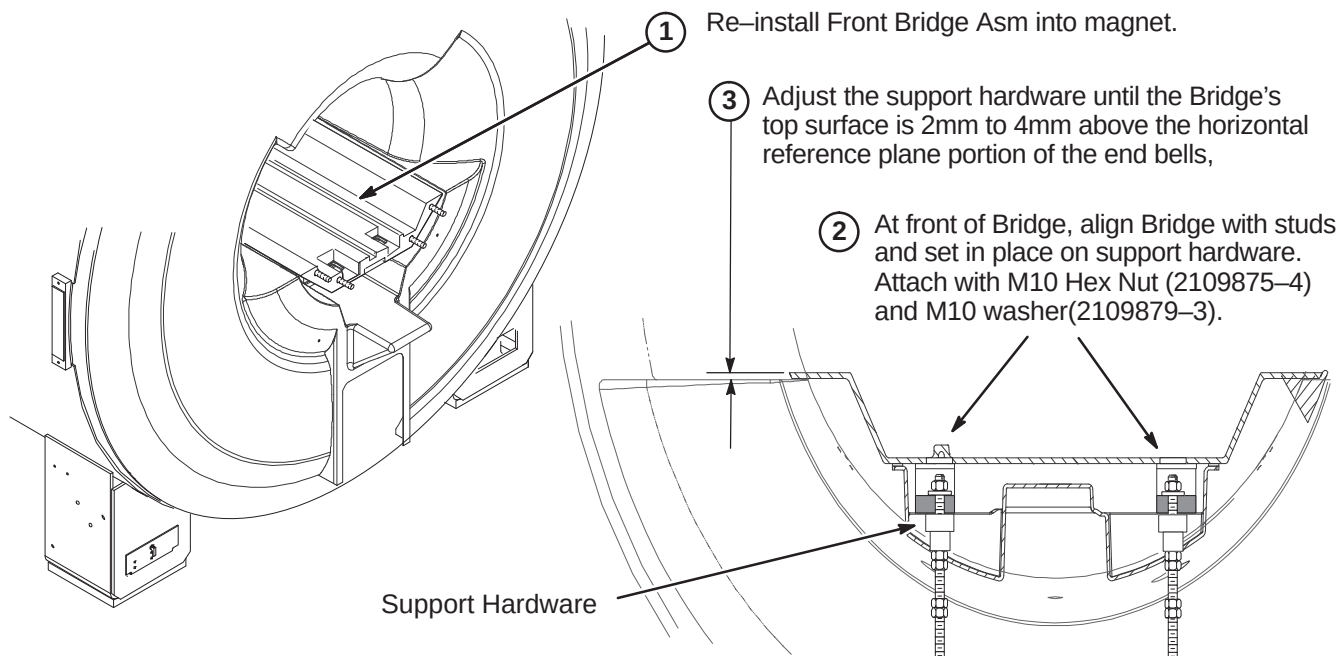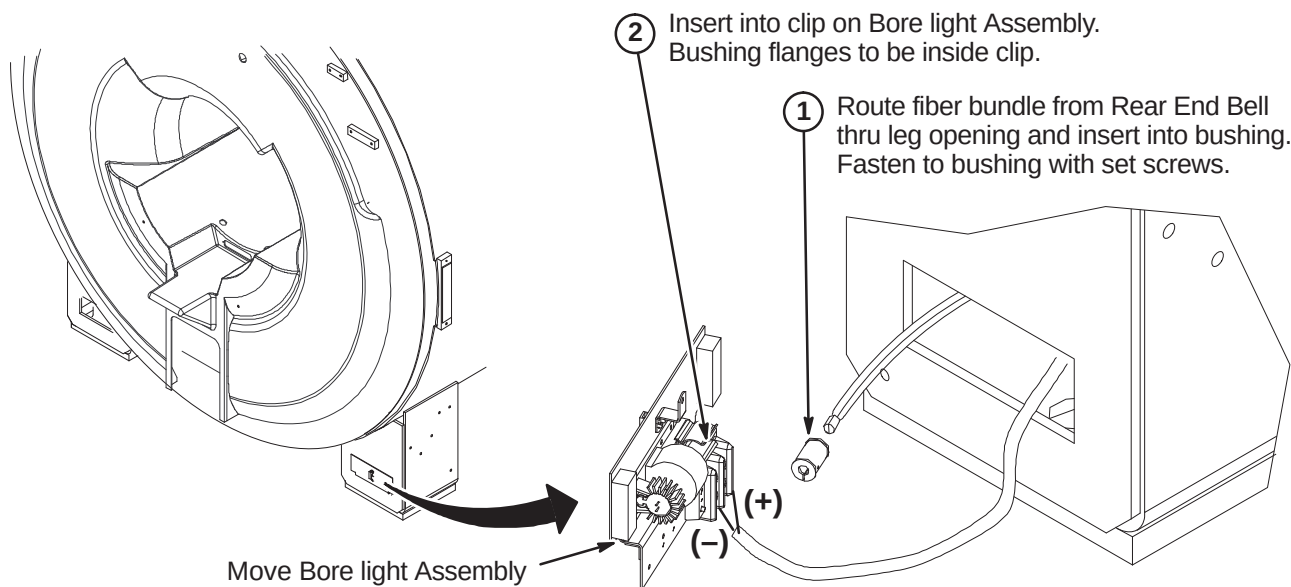4-5-20 Install Front Bridge**4-5-21 Connect Fiber Optic Bundle to Bore Lamp Assembly**

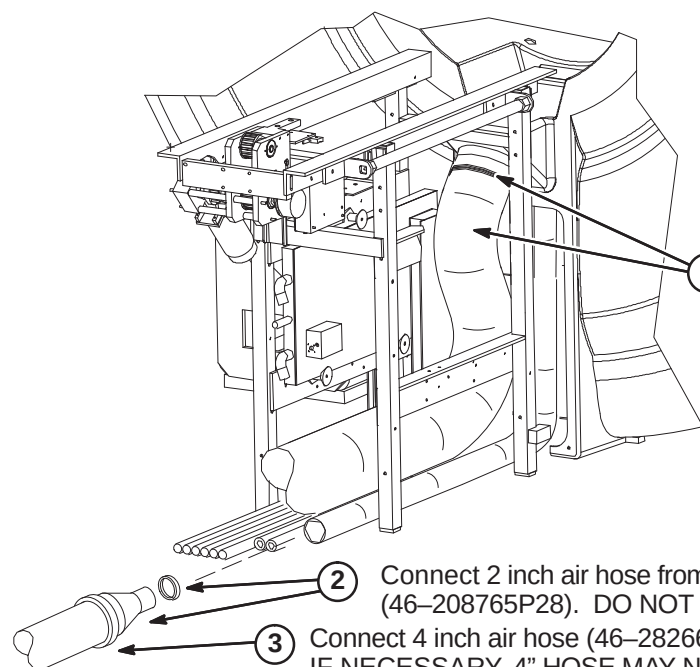4-5-22 Connect Air Hoses to Rear End Bell**Note:**

The following items are part of Blower Hose Kit, 2129412-6.

1 Connect 6 inch air hose (46-282666P14) to Rear End Bell with 6 inch clamp (46-208765P29).

2 Connect 2 inch air hose from Rear End Bell to Reducer (2175143) with 2 inch clamp (46-208765P28). DO NOT CUT OR REDUCE THE LENGTH OF THE 2" HOSE.

3 Connect 4 inch air hose (46-282666P13) to reducer with 4 inch clamp (46-208765P34). IF NECESSARY, 4" HOSE MAY NEED TO BE CUT TO REDUCE LENGTH.

4-5-23 200 Gauss Zone Limit for Movement of Blower box

Refer to the illustration and table below for applicable 200 Gauss zone for magnet at site.



THE BLOWER BOX IS HIGHLY FERROUS! IF MAGNET IS RAMPED, ANY MOVEMENT OF THE BLOWER BOX REQUIRES TWO PEOPLE AND MUST BE OUTSIDE THE 200 GAUSS ZONE.

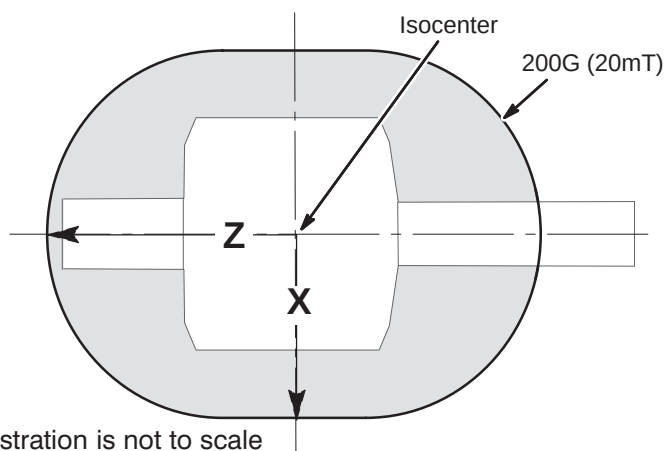


Illustration is not to scale

Approximate Distance from isocenter to 200G Line		
	X	Z
1.5T	150 cm (59 in) ¹	205 cm (81 in) ¹
Unshielded Magnet	Measure ²	
Note 1: These distances represent maximum distances for the magnet strength listed. Use gauss meter if needed. Note 2: For unshielded magnets, the magnetic field must be measured using a gauss meter.		


4-5-24 Installation of Blower Box and Air Hose from Rear of Magnet to Blower Box

CAUTION

***THE BLOWER BOX CONTAINS FERROUS MATERIAL. The proper mounting of the Blower Box is critical to safety. It MUST be mounted per these instructions.**

WARNING!

RF SHIELD INTEGRITY MUST BE MAINTAINED FOR MOUNTING BLOWER BOX WITHIN THE MAGNET ROOM

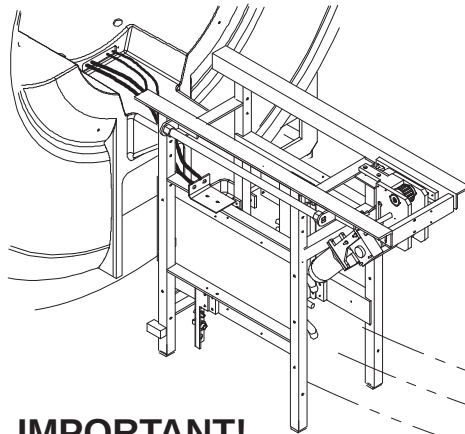**IMPORTANT!**


Location of cabinet and type of fastener used for securing Blower Box outside of minimum service area should be specified in the Site Design (architectural)/Construction drawings. If proper mounting area is not specified, contact the Project Manager of Installation (PMI).

- ⑥ Route 4 inch and 6 inch air hoses to Blower Cabinet. Trim to appropriate length and attach to vent openings with:

(1) 46-208765P34 4 inch Hose Clamp
 (1) 46-208765P29 6 inch Hose Clamp


- ① Refer to Site Design (architectural)/Construction drawings for locating and securing Blower Cabinet (2121199).
- ② Check for correct installation of the anchors:
 - Confirm that the anchors are physically secured to the building structure, not to a removable computer floor, raised floor, or false floor.*
 - Confirm there are at least four (4) anchors.
 - If any or the anchors appears wobbly or loose, stop the installation and contact the PMI to fix anchoring.
- ③ Remove any mounting hardware from the anchors.
- ④ Position the blower box on top of anchors.
- ⑤ Secure the blower box using the hardware provided by the RF shield room vendor to attach.

NOTE: If needed, Blower Hose Cover can be cut and modified.

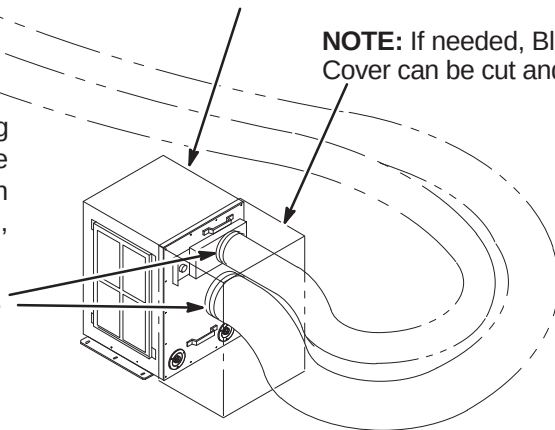

4-6 INSTALL REAR PEDESTAL

4-6-1 Install Isolation Filters to Hybrid Splitter

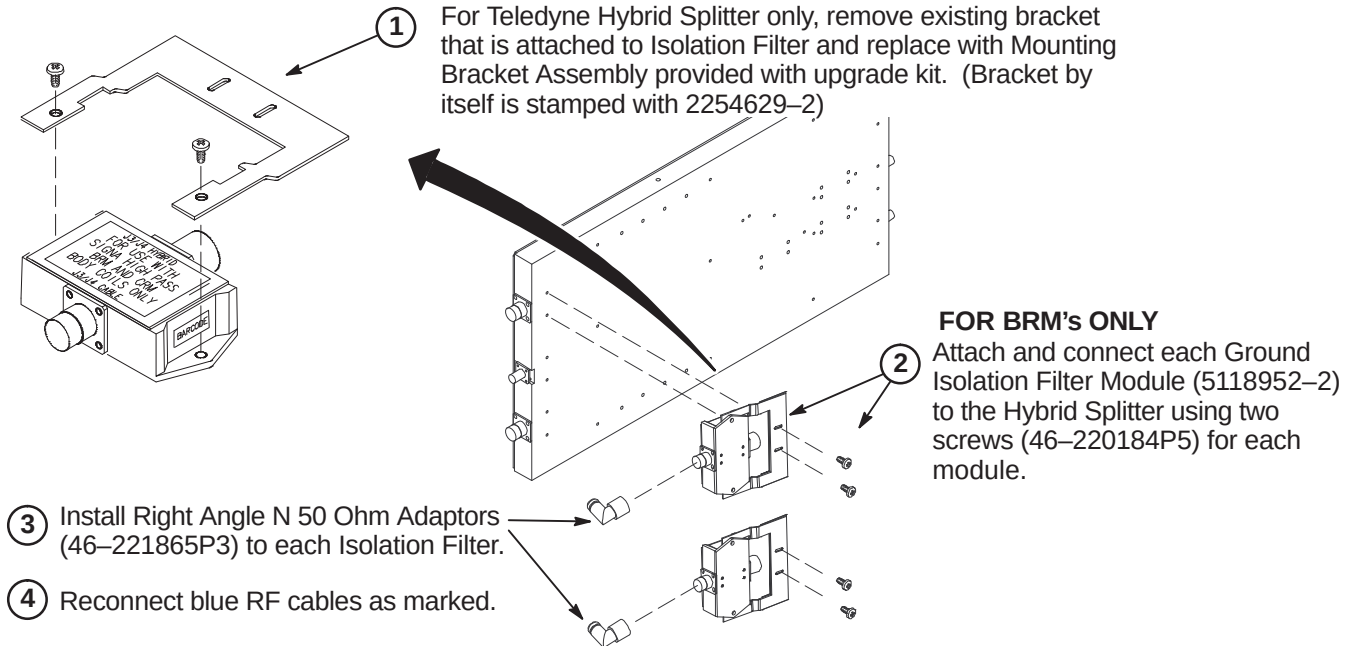


4-6-2 Connect Rear Pedestal to Magnet Enclosure

WARNING!

FERROUS MATERIAL HAZARD!


Rear Pedestal requires at least three people to remove from the magnet room. The longitudinal drive motor is ferrous. Keep the rear pedestal, especially the longitudinal drive motor, as far from the magnet as possible at all times.

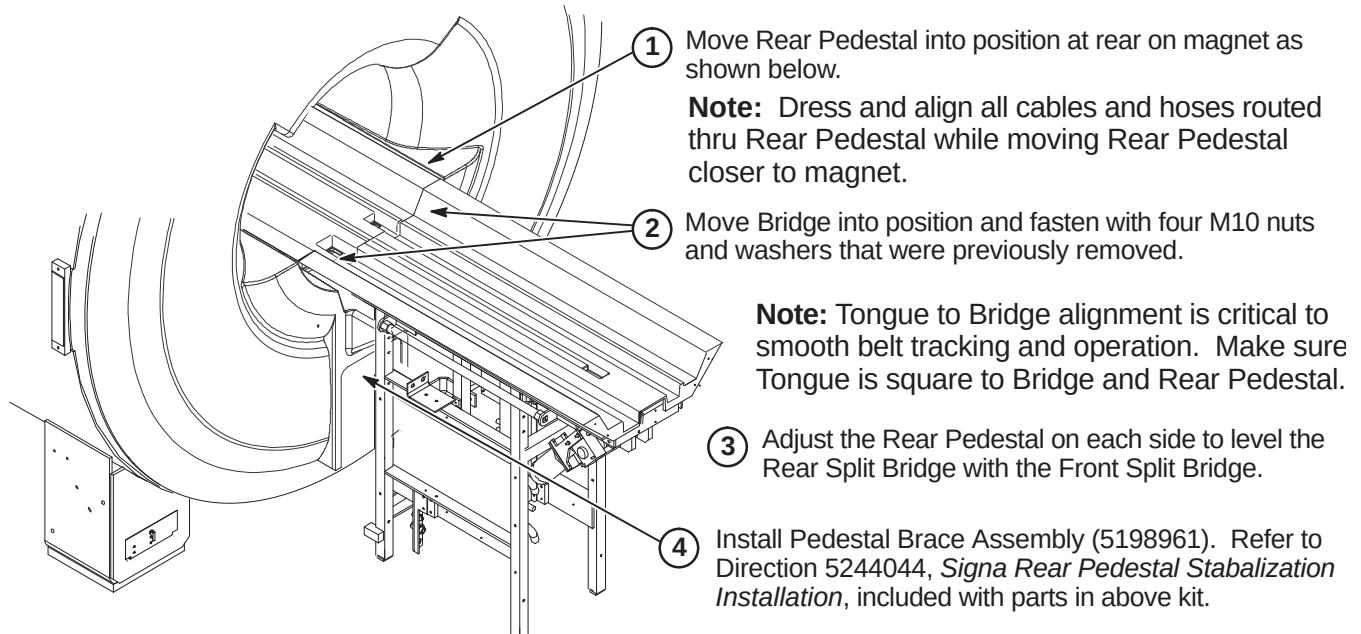



4-6-3 Leveling of Bridge

NOTE: Make sure that surfaces between front and rear Bridge surfaces are even. This will provide smooth travel for the Carriage wheels.

NOTE: Rear Pedestal must be adjusted to maintain full contact between bottom of Bridge and top of Rear Pedestal.

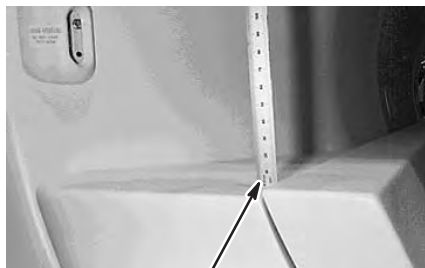

2mm to 4mm above End Bell surface.

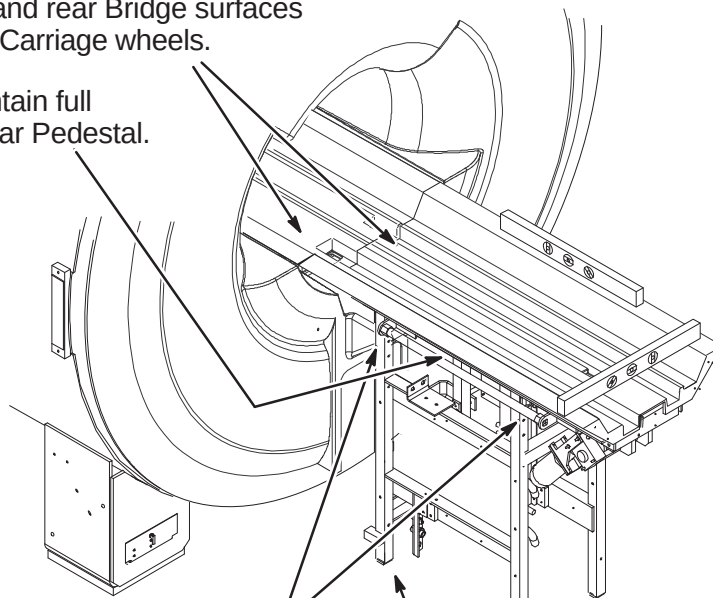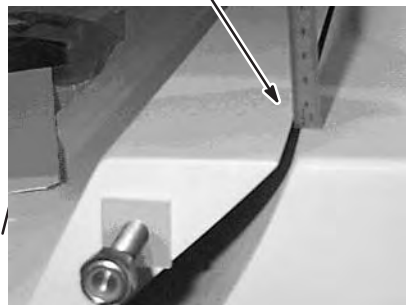


① Level Bridge and make sure that Bridge surface is no less than 2mm and no more than 4mm above the surface of End Bells, as indicated in photos at left. Tighten Vertical Adjustment Fasteners.

② Adjust leg screws on both sides as needed to level Rear Pedestal.

4-6-4 Height Adjustment At Front Of Bridge (If Needed)



2mm to 4mm above End Bell surface.

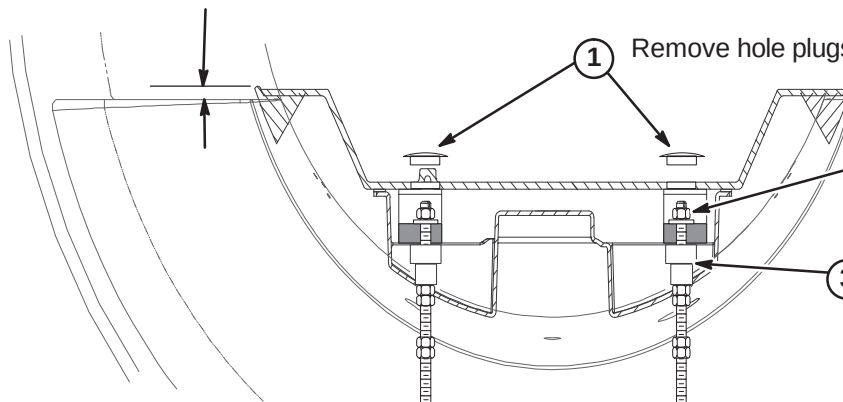

① Remove hole plugs from each side of Bridge.

② Remove M10 nut and washer from each stud.

③ Adjust Support Nut (2293381) as needed to achieve correct height.

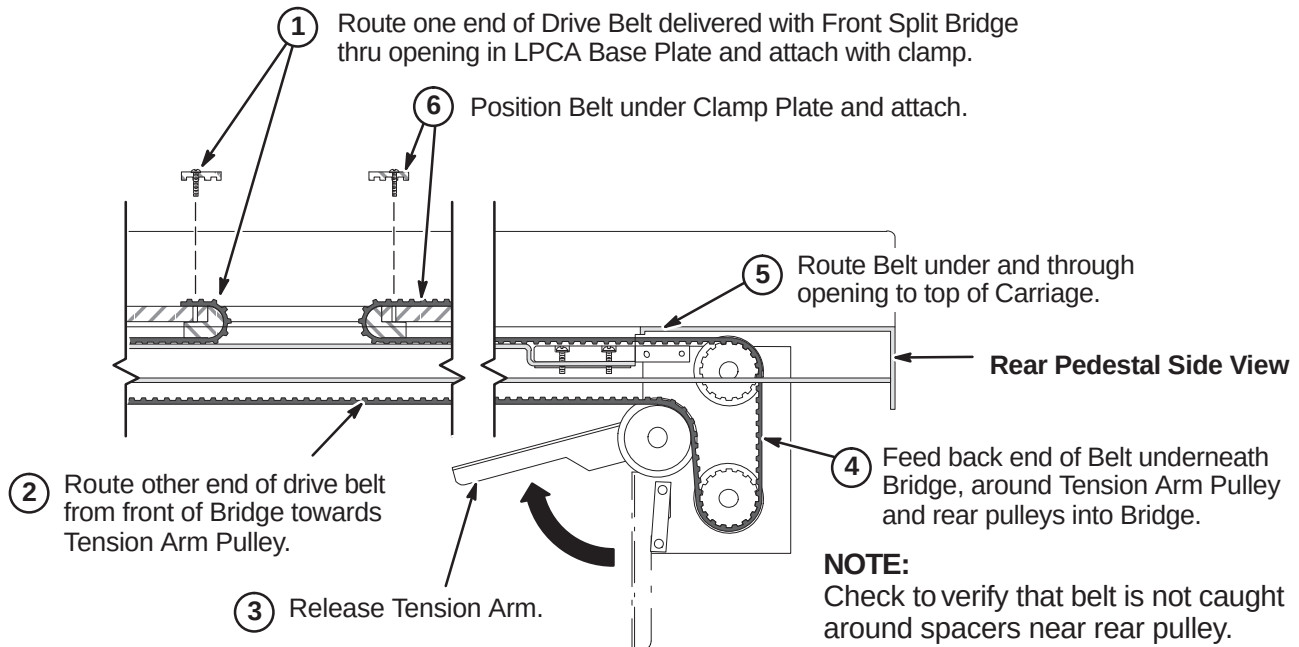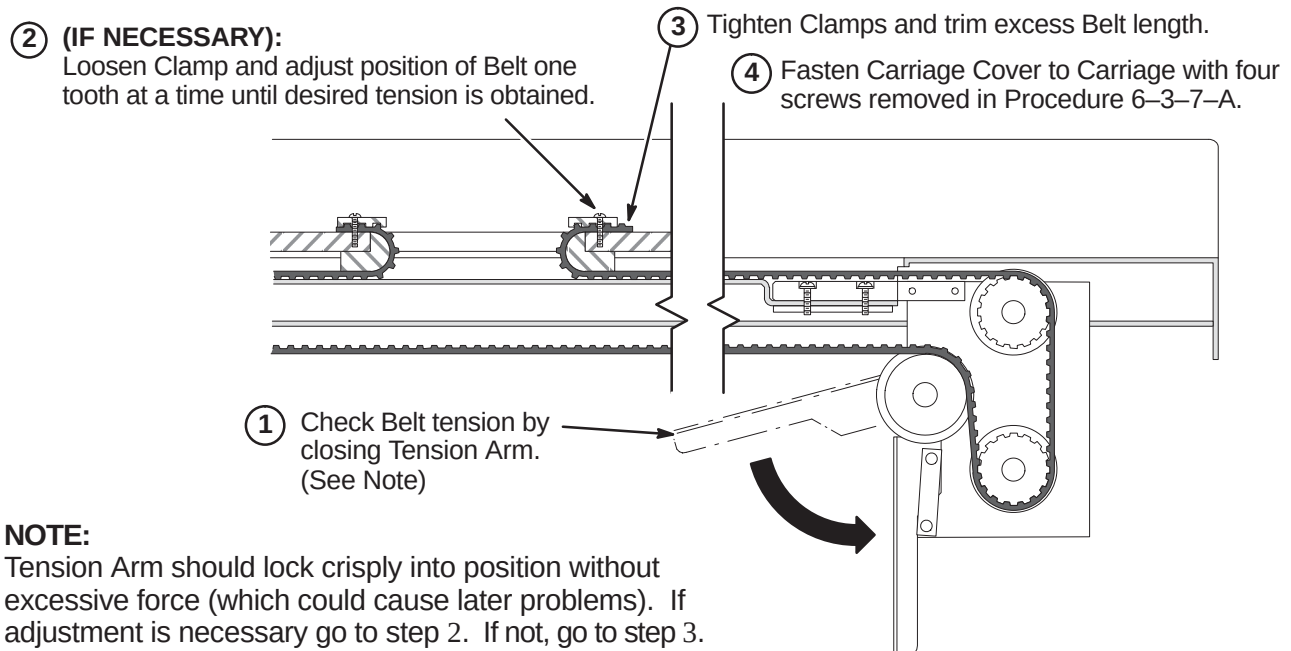4-6-5 Route and Attach Drive Belt**4-6-6 Final Adjustment of Drive Belt**

4-6-7 Installation of Sensor Flag


NOTE:

Function of the Flag is to interrupt the reflected limit switch sensor beam when the Carriage is in home position all the way forward in the bore.

NOTE:


The following Step is needed if flag was removed during routing of drive belt.

- ① Install Flag Assembly on Drive Belt in line with Bridge mounted Photo Sensor.

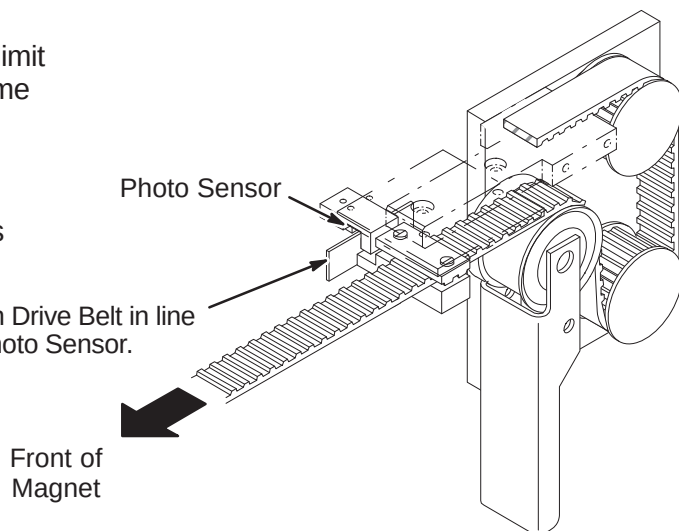


4-6-8 Installation of LPCA Cover

- ③ Attach the cover to the base plate using four screws previously removed.

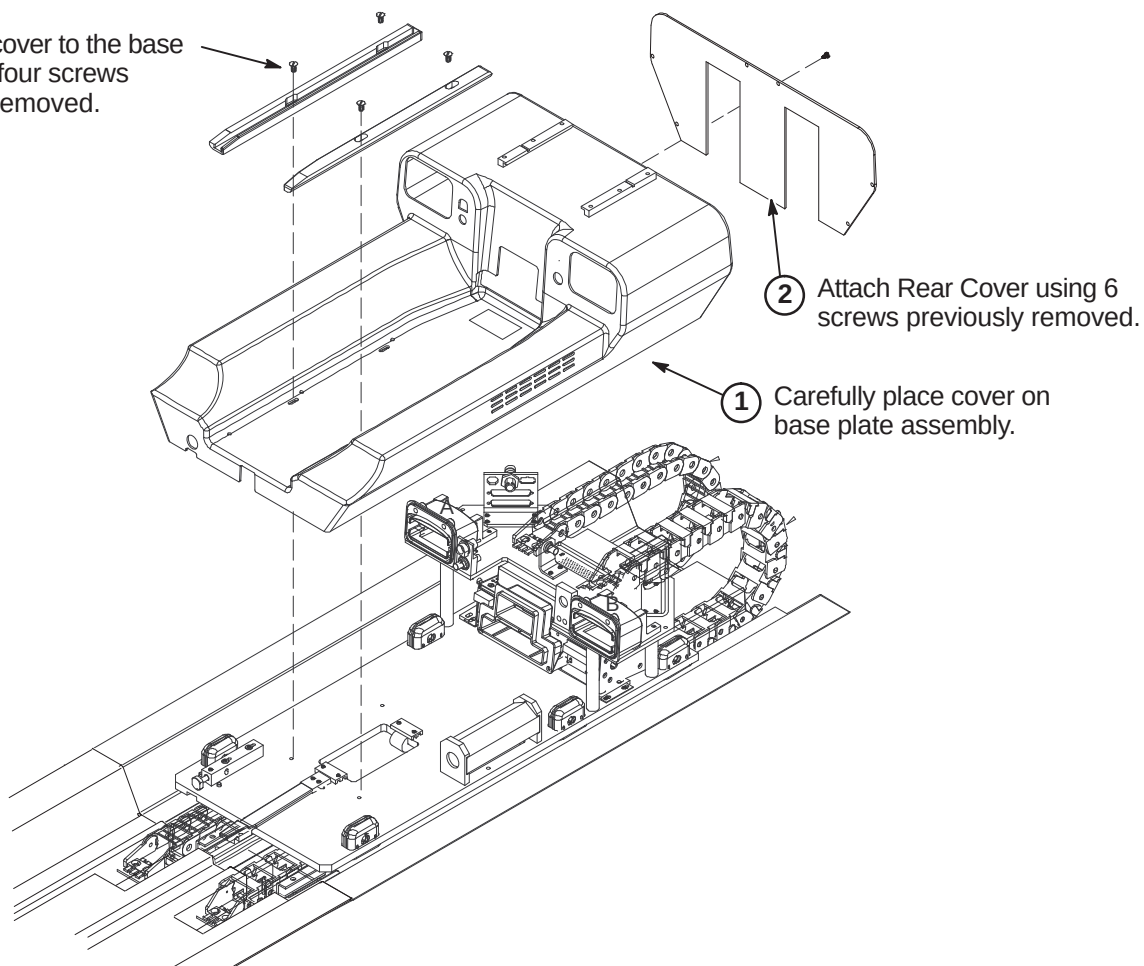


4-6-9 Reconnect Cables at Rear Pedestal

1. Route and connect RF Bias cables from RF Coil to bracket as marked.
2. Route and connect Run 781 from Bias cable J5 to J75 on the Penetration Panel.
3. Route and connect Run 782 from Bias cable J6 to J76 on the Penetration Panel.
4. Wrap both Bias connections with the previously removed 177 x 200 Connector Wrap.
5. Reconnect cables to Rear Pedestal components that were previously disconnected.
6. Route and connect I and Q RF Drive input cables from RF Coil and connect as marked to J1 and J2 on bracket.

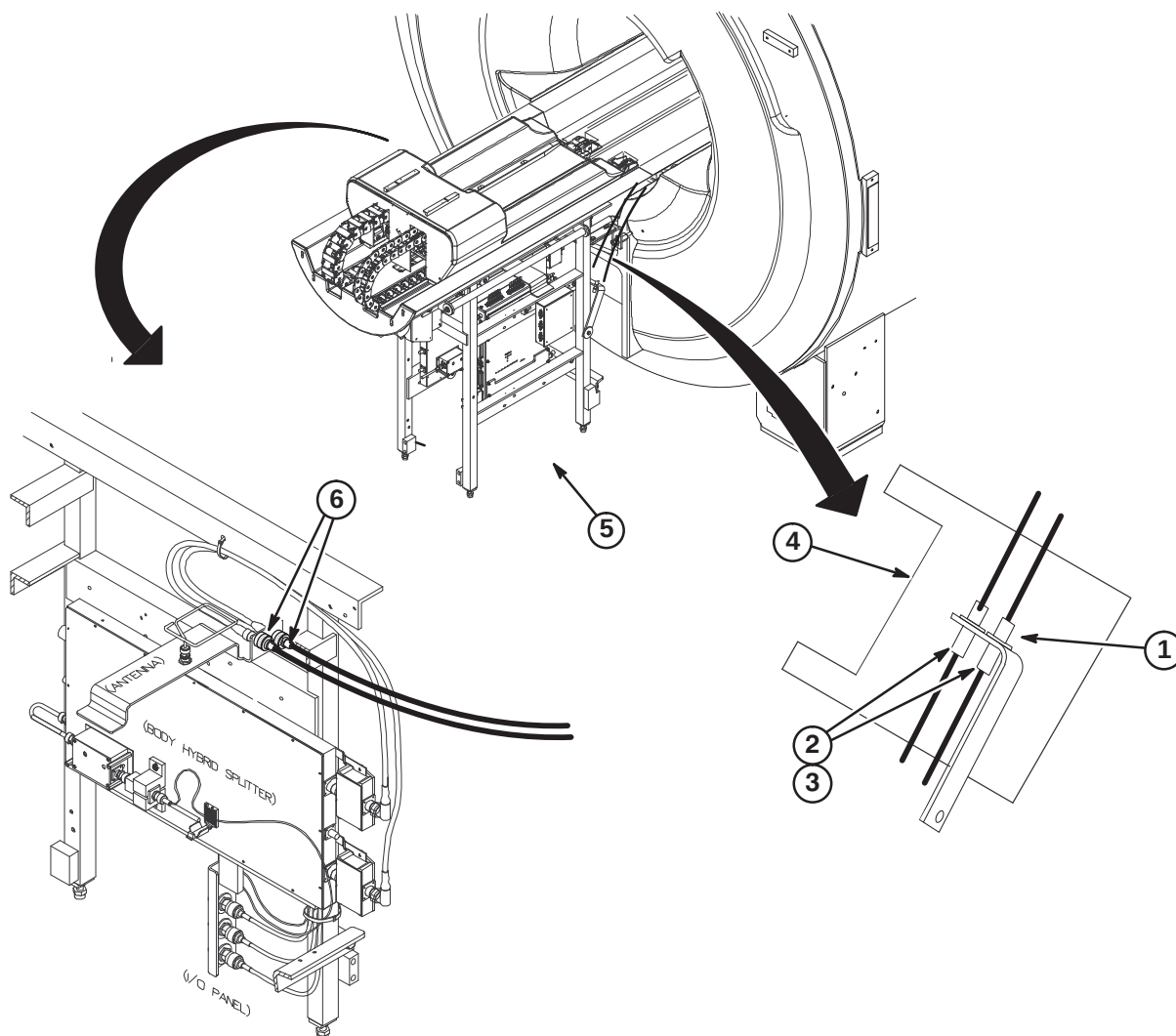

4-7 INSTALL CHILLER

4-7-1 Position Water Chiller Cabinet And Connect Hoses

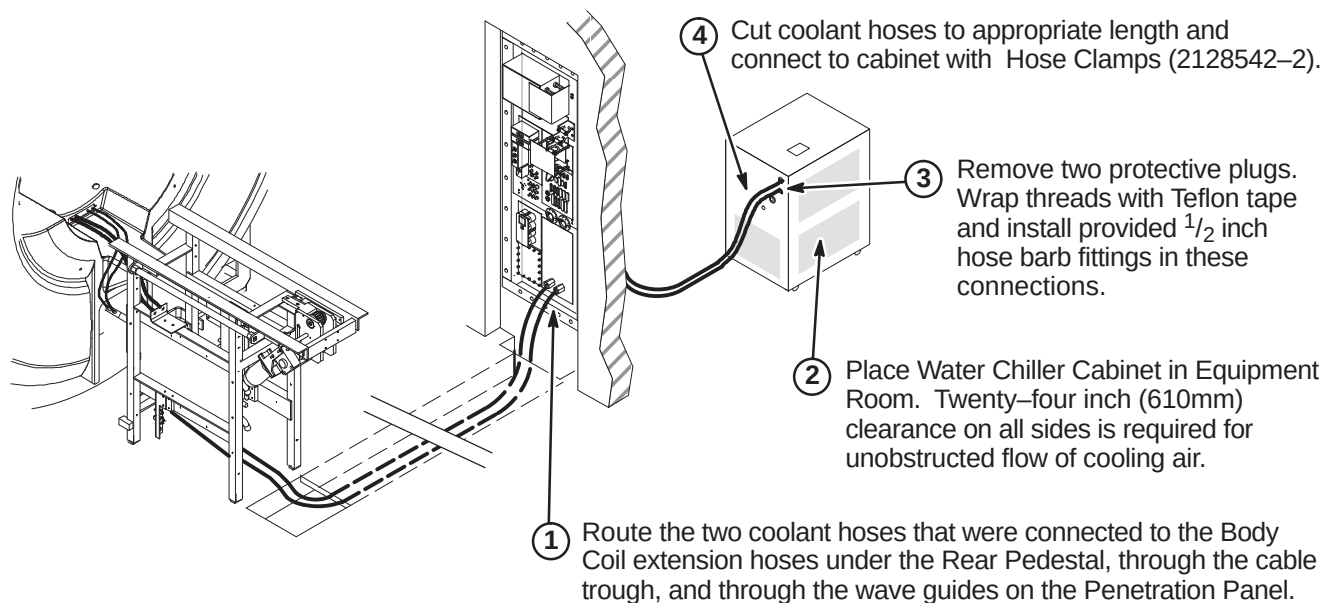



4-7-2 Install Power Cord

CAUTION

Do not start Water Chiller pump before filling reservoir with coolant.
Consult vendor Water Chiller manual –
Installation, Operation, and Service Manual that is shipped with the Water Chiller.

- ① Route and connect power cable Run 970 (2371516) to PDU as marked (TS1-16,17,GND).

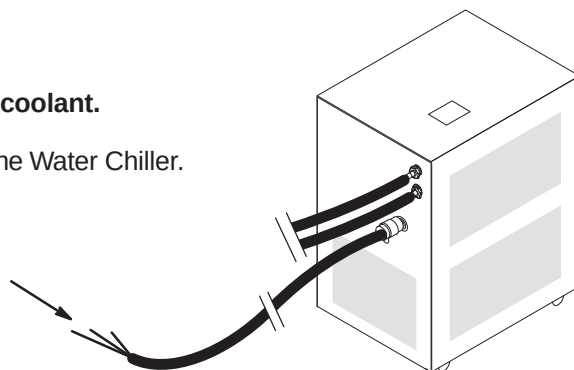

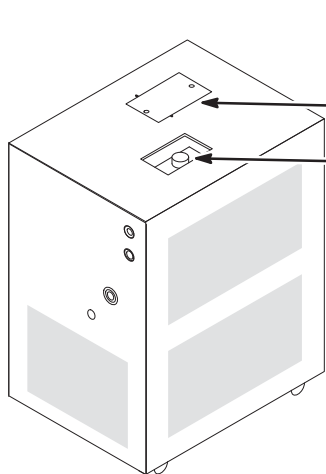4-7-3 Filling System With Coolant And Start Up

- ① Locate the twelve (12) gallons of Body Coil Coolant (2138791) that was shipped with System.
- ② Remove cover.
- ③ Consult vendor Water Chiller manual – *Installation: Filling Requirements*. Initial Filling of Reservoir, and the following steps for the proper procedure on priming and filling the system with coolant.
- ④ Consult vendor Water Chiller manual – *Operation: Start Up*. Follow instructions in manual for starting water chiller operation.
- ⑤ **Verify unit is set to control at 20° C.**
- ⑥ If supplied, attach plastic sleeve to side of Water Chiller and insert vendor manual. This will allow easy access to any future service/ maintenance procedures. If plastic sleeve is not present, attach vendor manual to Water Chiller.

4-8 INSTALL GRADIENT CABLES

Three cables have been delivered to site. All have three lead wires with terminals at each end. The leads on this cable end are connected to the Equipment room and Magnet room sides of the Gradient Filter on the Penetration Panel. This procedure cuts the cables to length and connects them to the HFD cabinet and the cable leads from the Gradient Coil.


Verify that 100 ft. (30.5 M) length of supplied Run 762, 763, and 764 cable is sufficient to route from HFD Cabinet to Penetration Panel and Penetration Panel to Body Coil Gradient Cables at rear of magnet. **Each cable has identical terminations at each end for connection to both sides of the Penetration Panel Gradient Filter.**

WARNING!


FERROUS MATERIAL HAZARD! THE RATCHETING CRIMP TOOL REQUIRED FOR THIS PROCEDURE CONTAINS FERROUS MATERIAL AND WILL BE STRONGLY ATTRACTED TO MAGNET AND MAY BECOME A DANGEROUS PROJECTILE. KEEP ALL FERROUS TOOLS AS FAR FROM THE MAGNET AS POSSIBLE.

IF MAGNET IS AT FULL FIELD – CUT GRADIENT CABLE TO LENGTH AND REMOVE FROM MAGNET ROOM TO PERFORM CABLE TERMINATION PROCEDURES.

4-8-1 Connect Gradient Cables to Equipment Room Side of Penetration Panel

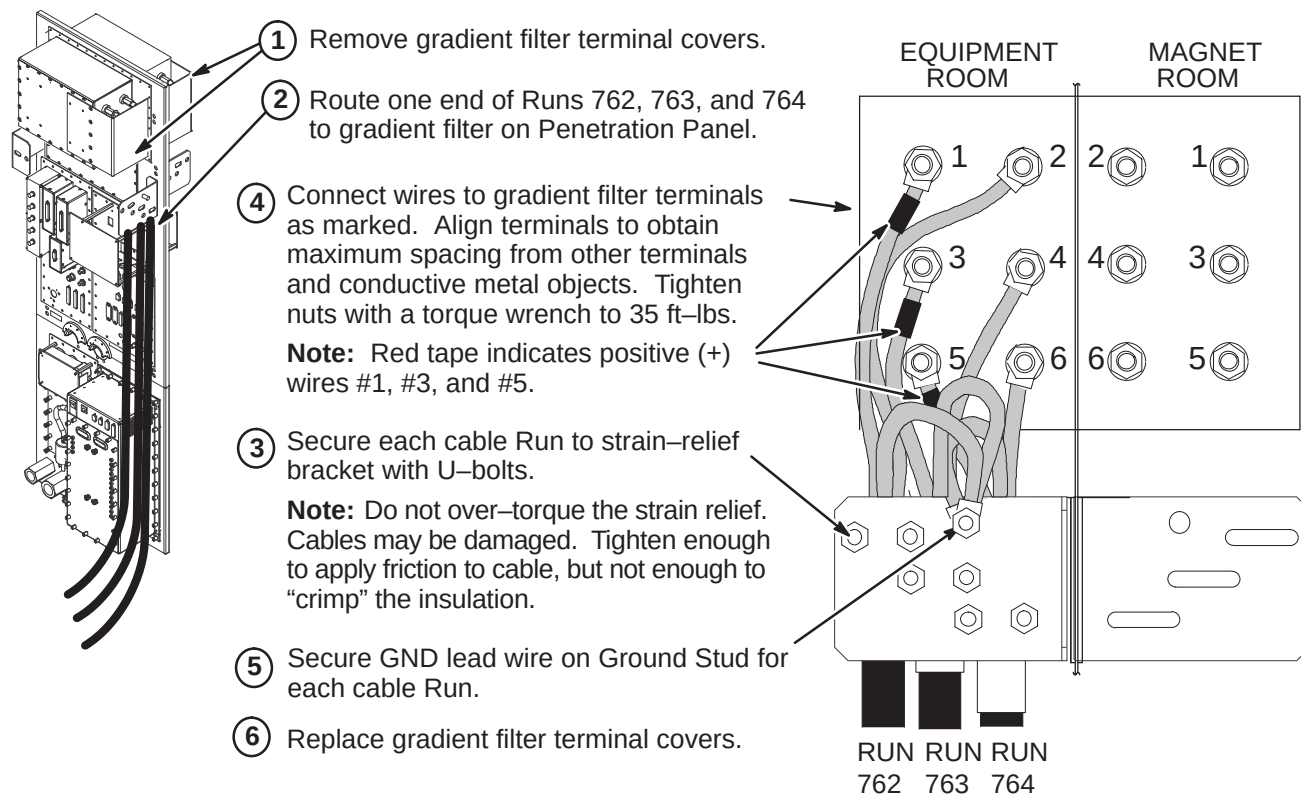

4-8-2 Route Runs 762, 763, and 764 to HFD Cabinet

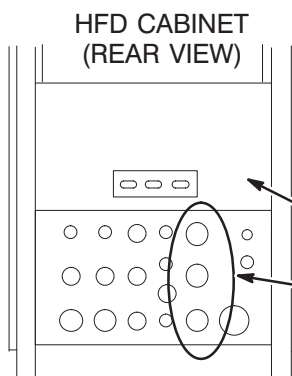


- ① Unroll and route each cable from equipment room side of Penetration Panel to HFD Cabinet. Verify sufficient slack of at least two feet (610mm) is provided for termination and connection at HFD Cabinet.
- ② Mark location for cutting.
- ③ Cut cables to length.
- ④ Remove Rear Cover Plate by removing seven screws that hold it in place.
- ⑤ Loosen screws on clamps for gradient output cables. After termination of cables in Procedure 4-8-3, insert ends of cables into Cabinet Interface as marked. (Runs 762 into J3, Run 763 into J4, and Run 764 into J5)

4-8-3 Terminate Runs 762, 763, and 764 for HFD Cabinet

- ① Refer to above table to determine length of jacket and foil to be removed from end of cable Runs 762 (X), 763 (Y), and 764 (Z).


	Run 762	Run 763	Run 764
Jacket Length to be Removed	8.0 IN (203mm)	11.0 IN (280mm)	14.0 IN (356mm)

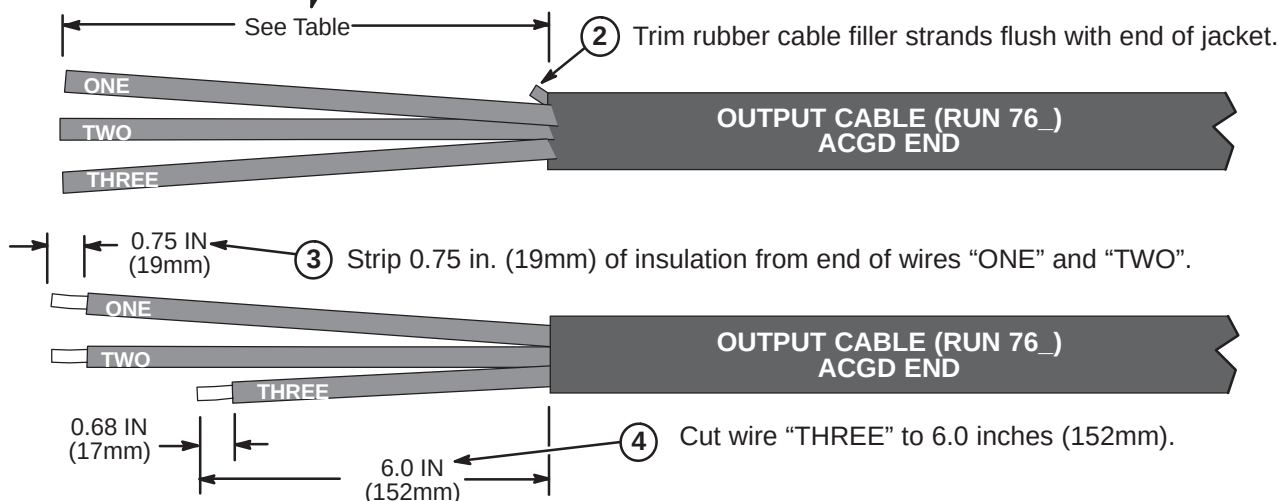

Note:

Steps below show completion of Run 762. Repeat same process for Runs 763 (Y+, Y-) and 764 (Z+, Z-).

Red tape, Y+ (RUN 763) and Z+ (RUN 764) labels are attached to wire "ONE".

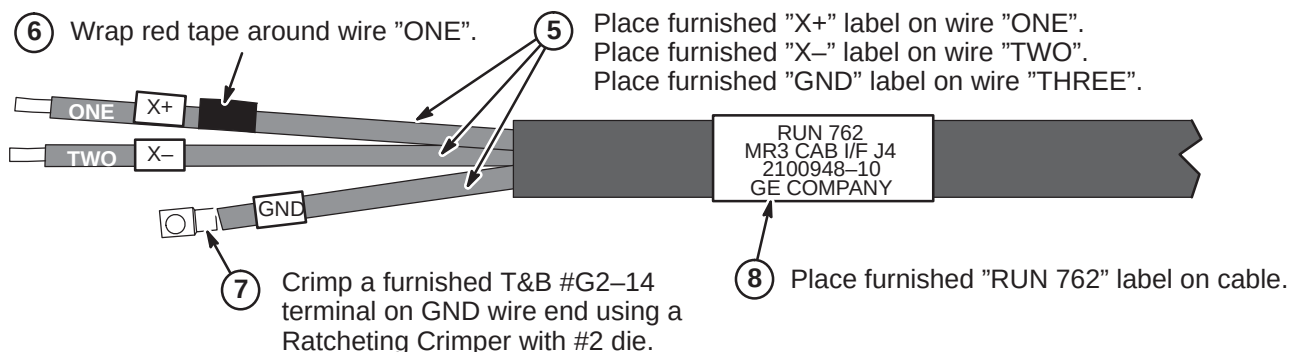

4-8-4 Connect Runs 762, 763, and 764 at HFD Cabinet


1. Connect gradient output cables with red-taped (+) #1 wires to + output terminals and black (-) #2 wires to - output terminals.

NOTE: MAKE SURE RUNS 762, 763, AND 764 DO NOT COVER 480V TERMINAL BLOCK


2. Tighten screws to 52–65 inch-lbs (5.88–7.34 N-M). This is about what can be applied using only your thumb and little finger on wrench. Loosen 1/2 turn and re-tighten to the same torque. Repeat approximately 6 times or until wrench returns to the same spot after torquing.

NOTE: OVER TIGHTENING OF SCREWS IN STEP 2 CAN DAMAGE TERMINAL BLOCK.


3. Route the GND leads of each run to the Ground Bus Bar studs. Connect from bottom to top.
4. Tighten clamps on Interface Panel to secure gradient output cables.

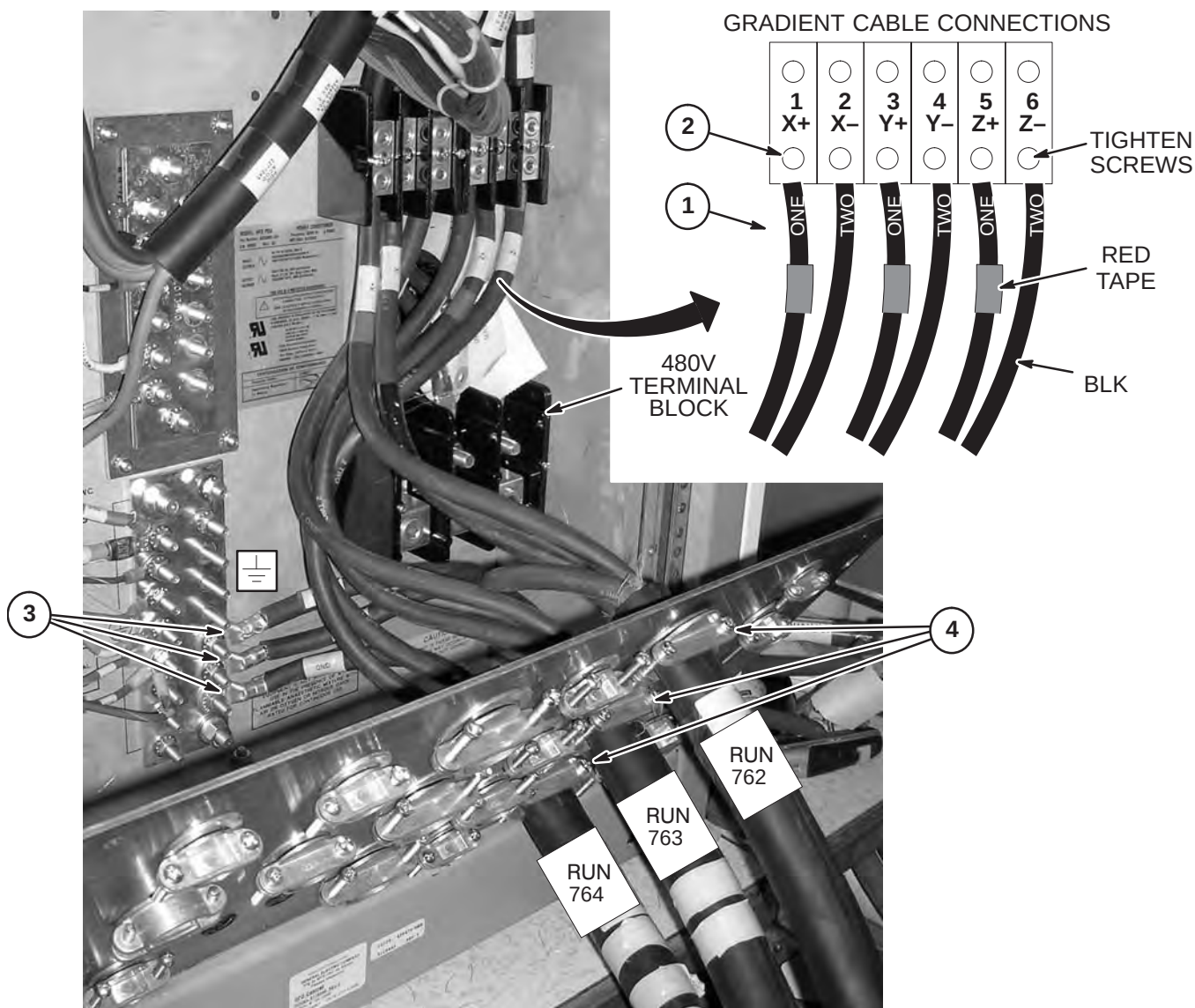


4-8-5 Connect Gradient Cables to Magnet Room Side of Penetration Panel

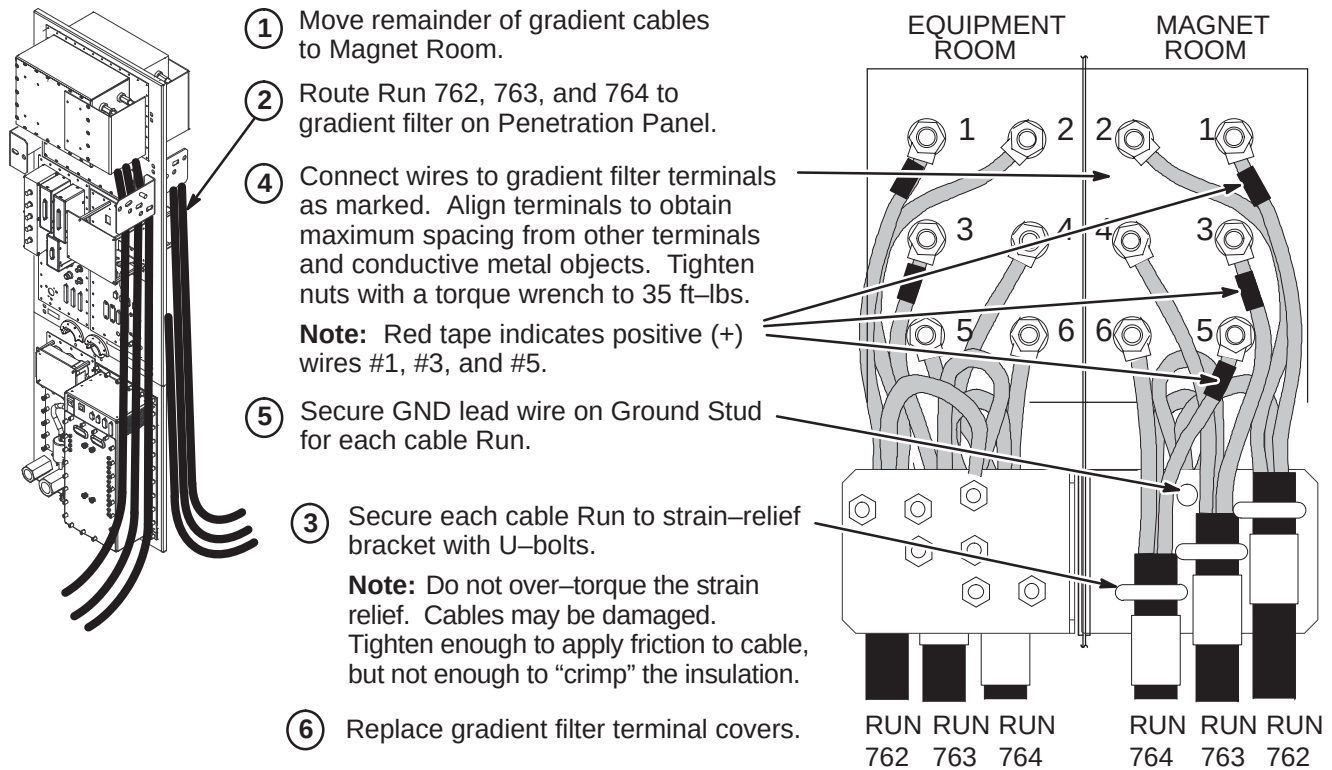

4-8-6 Route Runs 762, 763, and 764 to Rear Pedestal Area

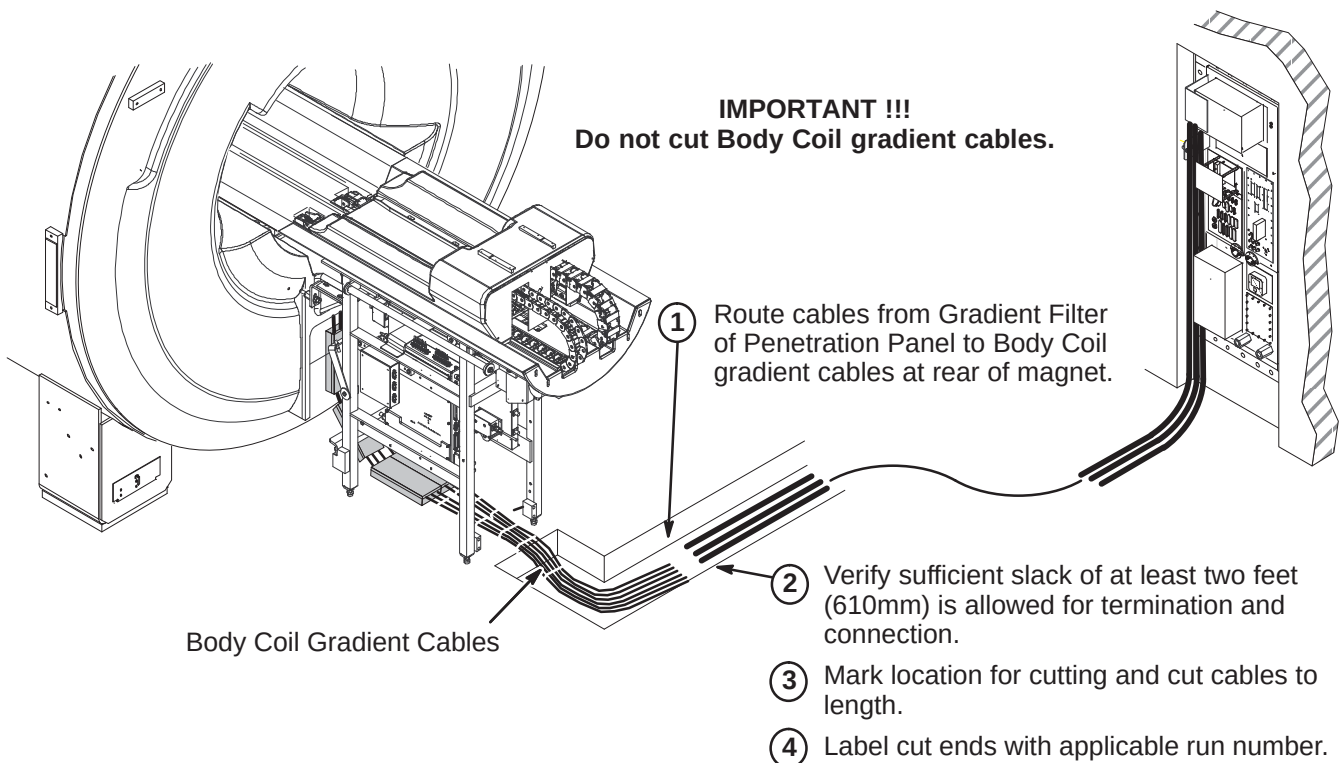

4-8-7 Terminate Runs 762, 763, and 764 for Magnet Room

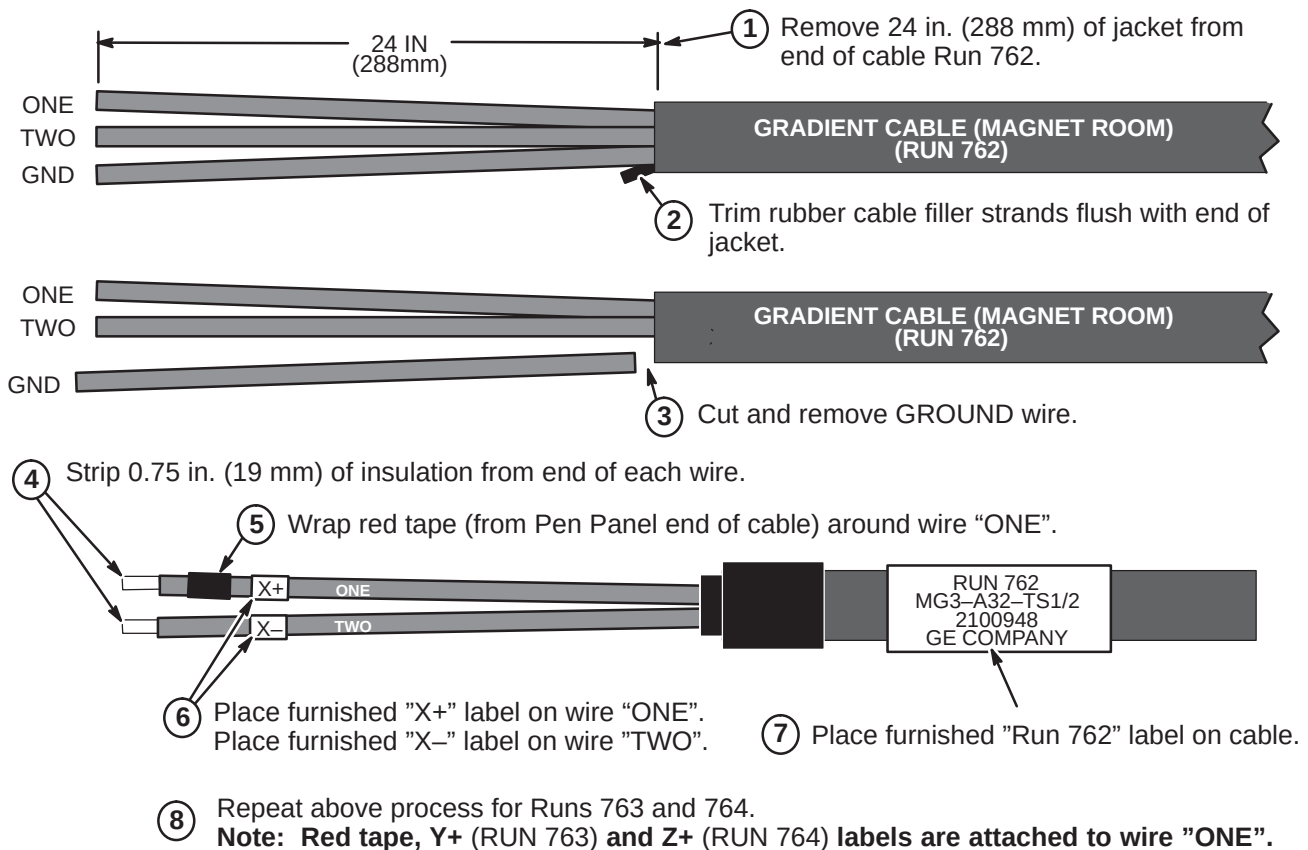


4-8-8 Connection of Body Coil Gradient Cables to runs 762, 763, and 764

WARNING!

FERROUS MATERIAL HAZARD! THE RATCHETING CRIMP TOOL REQUIRED FOR THIS PROCEDURE CONTAINS FERROUS MATERIAL AND WILL BE STRONGLY ATTRACTED TO MAGNET AND MAY BECOME A DANGEROUS PROJECTILE. KEEP ALL FERROUS TOOLS AS FAR FROM THE MAGNET AS POSSIBLE.

IF MAGNET IS NON-ACTIVE SHIELD AND IS AT FULL FIELD – MOVE END OF BODY COIL EXTENSION CABLES TO FARTHEST POINT AWAY FROM MAGNET TO PERFORM BUTT SPLICE CONNECTION TO RUNS 762, 763, AND 764 COMING FROM PENETRATION PANEL.

4-8-8 Connection of Body Coil Gradient Cables to runs 762, 763, and 764 (Continued)

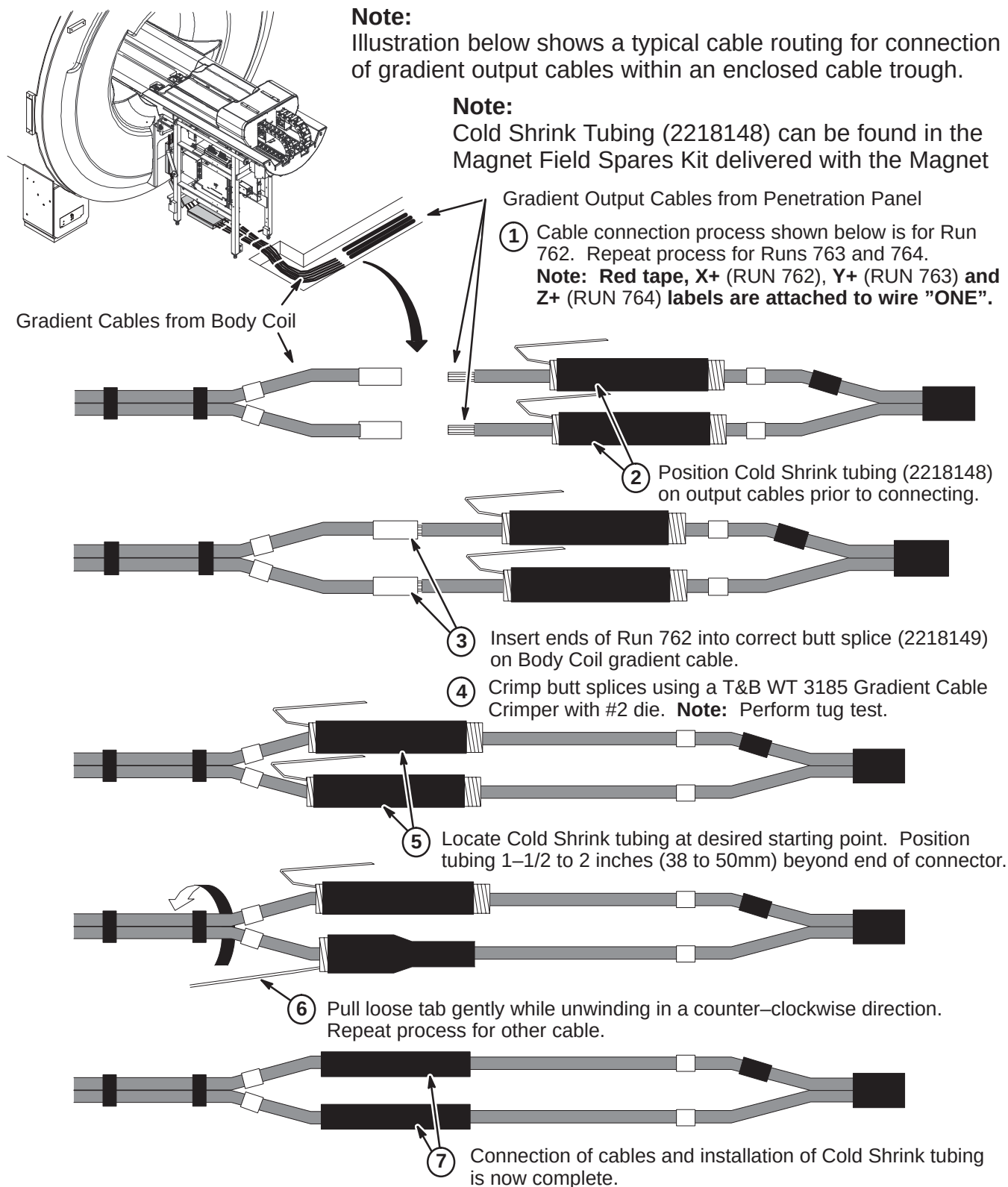

4-9 FINAL ENCLOSURE INSTALLATION


4-9-1 Install Enclosure Top Covers

Note: All enclosure seams and gaps visible to the operator and patient shall be less than 0.25 inches (6.35mm). This applies to the mating seam of the two Side Top Covers to the Front Top Cover


Gap uniformity along any continuous seam shall be ± 0.0625 inch (1.59mm)

- 1 Install front and rear Top Covers as shown using M10 x 25mm Low Profile hex screws and blue Loctite.

- 2 Place a level against vertical surface of Top Covers. If top of covers are leaning towards magnet, adjust Cosmetic Cover Ring and/or left and right Arc Covers.

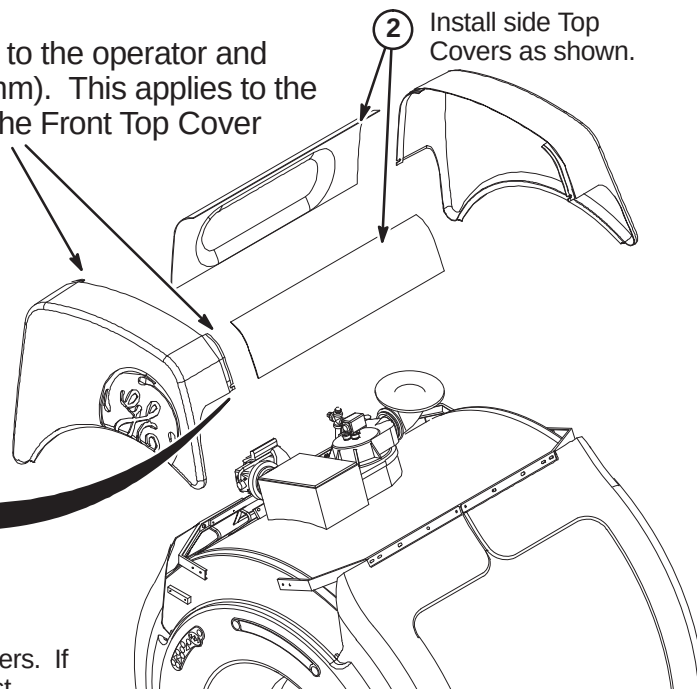

4-9-2 Install Base Trim Covers

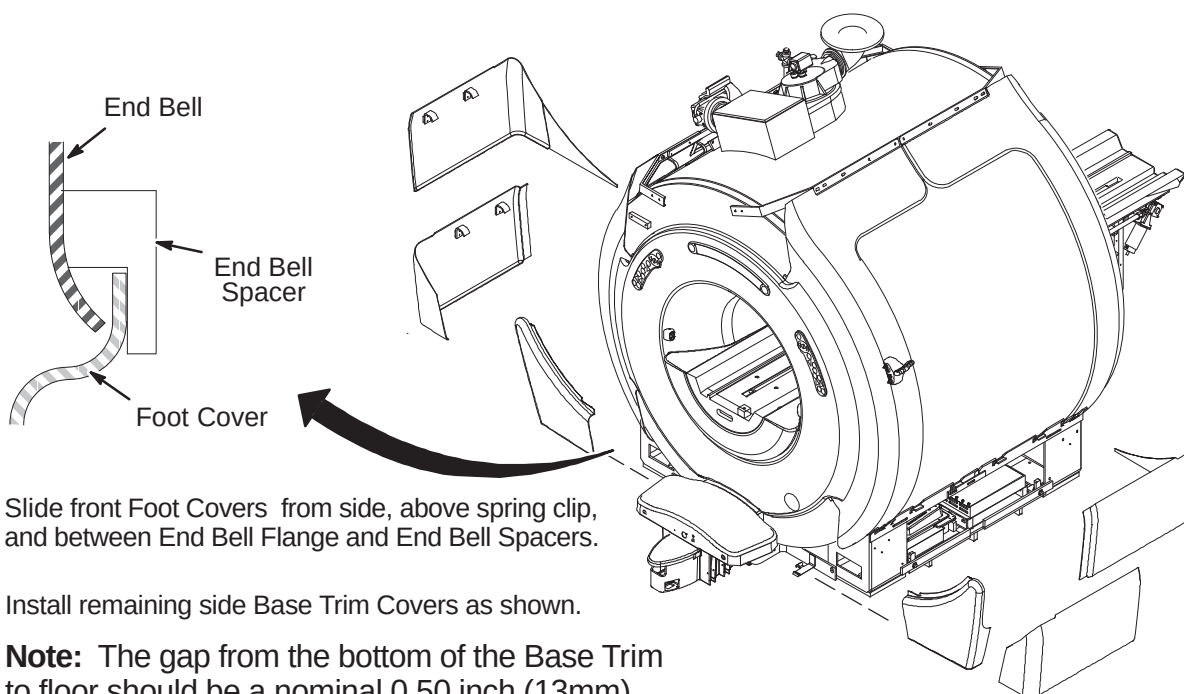


- 1 Slide front Foot Covers from side, above spring clip, and between End Bell Flange and End Bell Spacers.
- 2 Install remaining side Base Trim Covers as shown.

Note: The gap from the bottom of the Base Trim to floor should be a nominal 0.50 inch (13mm).


4-9-3 Install Rear Pedestal Covers

- ② Pedestal Covers (2217704 & 2217705) require an opening to be cut out depending on routing of cables and hoses. Covers are marked in three places (left side, right side, or rear) for cutting.
- ③ After cutting out appropriate opening, hang covers on Rear Pedestal.
- ④ If covers need adjustment, loosen the screws at each end of the support bar and reposition bar as necessary.

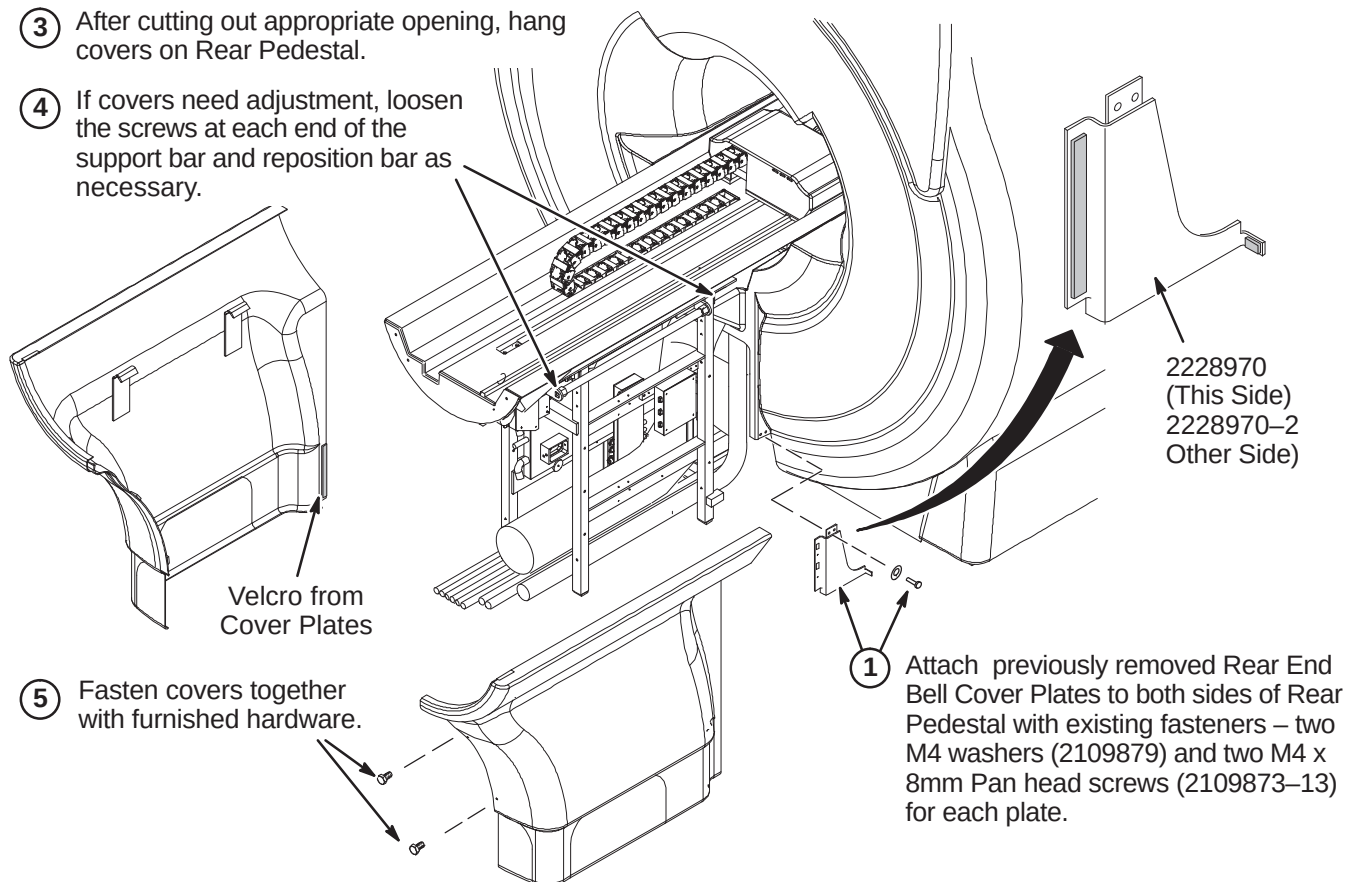

4-9-4 Install Rating Plates


- ① If not already applied, place the following Rating Plates on Bridge Rear Plate:
 46-302200P12 Kit Added Rating Plate – No S/N
 2142857 ETL List for Global Certification
 2117390-2 CE International Approval
 2132620 CE Marking Label Class A Group 1

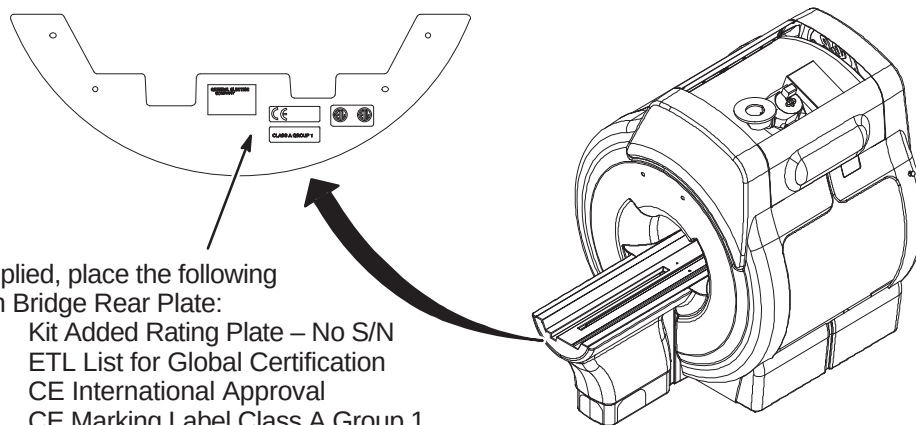

4-9-5 Install Laser Alignment Light Warning Labels

A sheet of Laser Alignment Light Warning Labels (2398124-2) has been supplied with the magnet, attached to the magnet bridge.

The magnet enclosure is shipped with a label in English attached below each of the laser lights.

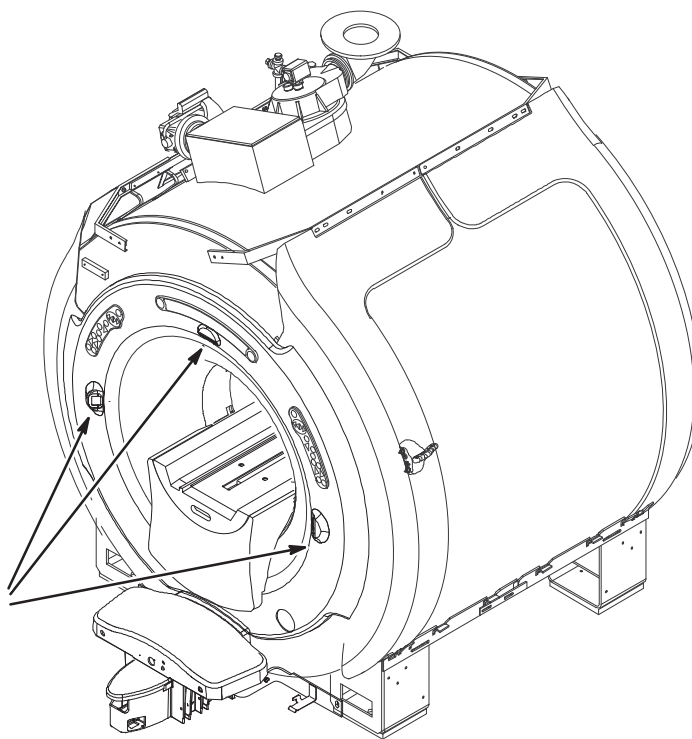

- ① For those sites which require a label other than English, peel off the appropriate label and affix over the English labels at the laser light locations.



French



German



Portuguese



Spanish




Italian



Swedish



Japanese

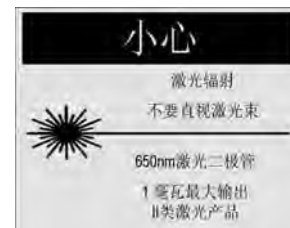

Chinese

This page left blank intentionally.

SECTION 5 – UPGRADE COMPLETION

5-1 MECHANICAL INSTALLATION COMPLETION

The following steps need to be reviewed and completed by GE Field Engineer, especially if the system was installed by non-GE personnel. Check off boxes at the front of each item below when complete.

- ☐ Check mechanical flowchart to make sure all steps are complete (Section 2)
- ☐ Resolve shipment shortages
- ☐ Resolve omissions by mechanical contractors
- ☐ Have the on-site GE Healthcare Field Engineer, Service Engineer, or Project Manager of Installation (PMI) to confirm that the blower box has been anchored properly and installed outside the minimum service area.
- ☐ Make sure all rating plates are installed
- ☐ Recheck all wiring connections using cable map (Section 4)
- ☐ Correct for wiring errors if necessary
- ☐ Complete Sections 5-2, 5-3, and 5-4.

5-2 RETURNING SHIPPING EQUIPMENT

The shipping material must be returned **as soon as possible**. Prompt return is required so that production shipments can continue with these labor saving dollies. These can be included with the obsolete equipment and materials being returned in Section 5-3. Recycling Operation will route them to appropriate location.

5-3 MATERIAL DISPOSITION GUIDANCE

Refer to Section 1-7. The equipment or material being disposed by the upgrade have a potential value for service of the installed base. The disposition of such must be controlled by GE. Disposition of the removed material is per GE Medical Systems Policy and Procedures.

5-4 FINAL COORDINATION

1. Replace cabinet covers that have not been previously replaced.
2. Replace Magnet Enclosure covers and Rear Pedestal covers (If required)
3. Dock Patient Transport in position at front of bore.
4. Record and enter applicable data into applicable site configuration files and records.
5. Complete Product Locator information. See Section 1-8.

Note

Failure to fill out and return Product Locator Cards may result in failure of your site to receive future FMI's.

6. Leave the site's set of Manuals and CD-ROM at site. Organize and set up reference cabinet for the Manuals.
7. Locate any new Material Safety Data Sheets (MSDS). They must be retained on site. Customer is to be informed that material with MSDS information was brought into site and customer should know/decide where MSDS be retained.
8. Verify that coordination of application support for instruction of site users on Software Release 14.x has been made before returning system to the users. Turn site over to applications who will instruct users.

5-5 POWER ON

1. Notify field service and other installation personnel that are working at the site that the PDU main disconnect locks and tags will be removed so PDU power can be turned on.
2. Restore power to the PDU from main disconnect.
3. Refer to the CD-ROM, Installation Calibration Flow, and perform the following:
 - PDU Voltage Checks Compact > PDU E-Stop Check
 - Startup Checks Compact
 - Transformer Taps Compact
4. Press the FULL ON button on the PDU front control panel.
5. Except for the HFD Gradient power modules, switch all cabinet circuit breakers to ON. The HFD Gradient power module circuit breakers should remain OFF for one hour after HVAC power up to allow room humidity to stabilize.
6. **Make sure the UPS circuit breaker is switched to the ON position.** This is required to feed power to the Operator Workspace.

SECTION 6 – SOFTWARE LOAD PROCEDURE

1. Power up system.
2. Log in as **root**, password is **operator**
3. Open Guided Install and change gradient coil to BRM
4. Set the ICW to “Install Mode” by
 - a. In C-shell type **su root**, password is **operator**
 - b. Type **cd /export/home/service/license**
 - c. Type **rm .READYTOGO**
 - d. Type **rm .inMaintenance** (This file may or may not be present)
 - e. Type **cd /export/home/service/log**
 - f. Type **rm serial.log**
5. Perform system calibrations per ICW
6. Perform SaveInfo.