



HIS/RIS Setup and Troubleshooting

[\(Click here for Summary Version\)](#)

1 Personnel Requirements

Required Persons	Preliminary Reqs	Procedure	Finalization
1	Not Applicable	15 mins	Not Applicable

2 Overview

This procedure documents the set up and calibration of the system to a clinical hospital information system (HIS) or a radiology information system (RIS) .

3 Preliminary Requirements

3.1 Required Conditions


Condition	Reference	Effectivity
Acquire HIS/RIS, filming network information from the customer administrator.	-	-

4 Procedure

4.1 Setup

1. From the Common Service Desktop, select the [Configure] button, select the [Guided Install] button, and in FE Mode login as *root* (password is *operator*).
2. Under [Connectivity], select [HIS/RIS].

Illustration 1: HIS/RIS Screen

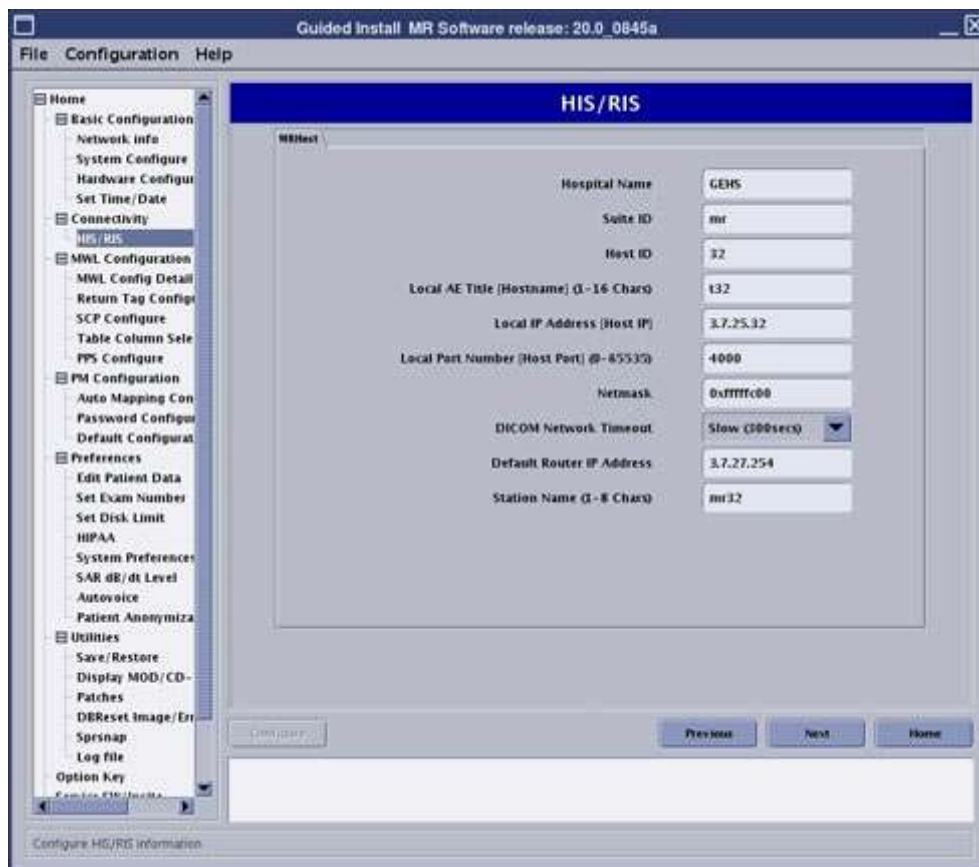

3. Enter applicable information in all fields (as provided by site) and click the [Configure] button.
4. Under [MWL Configuration], select [MWL Config Details].

Illustration 2: MWL Config Details Screen

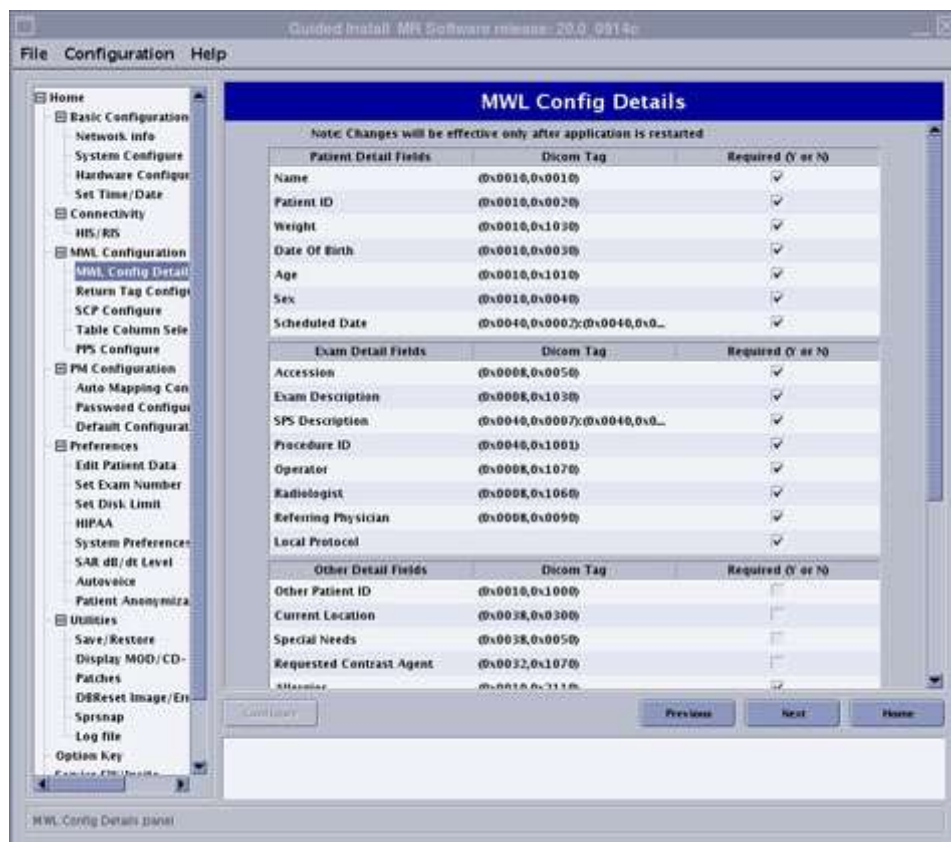



Illustration 3: MWL Config Details Screen (Continued)



5. If any of fields, other than the defaults, are desired to appear in the Modality Worklist, select or deselect appropriately and click the [Configure] button.

6. Under [MWL Configuration], select [Return Tag Configure].

Illustration 4: Return Tag Configuration Screen

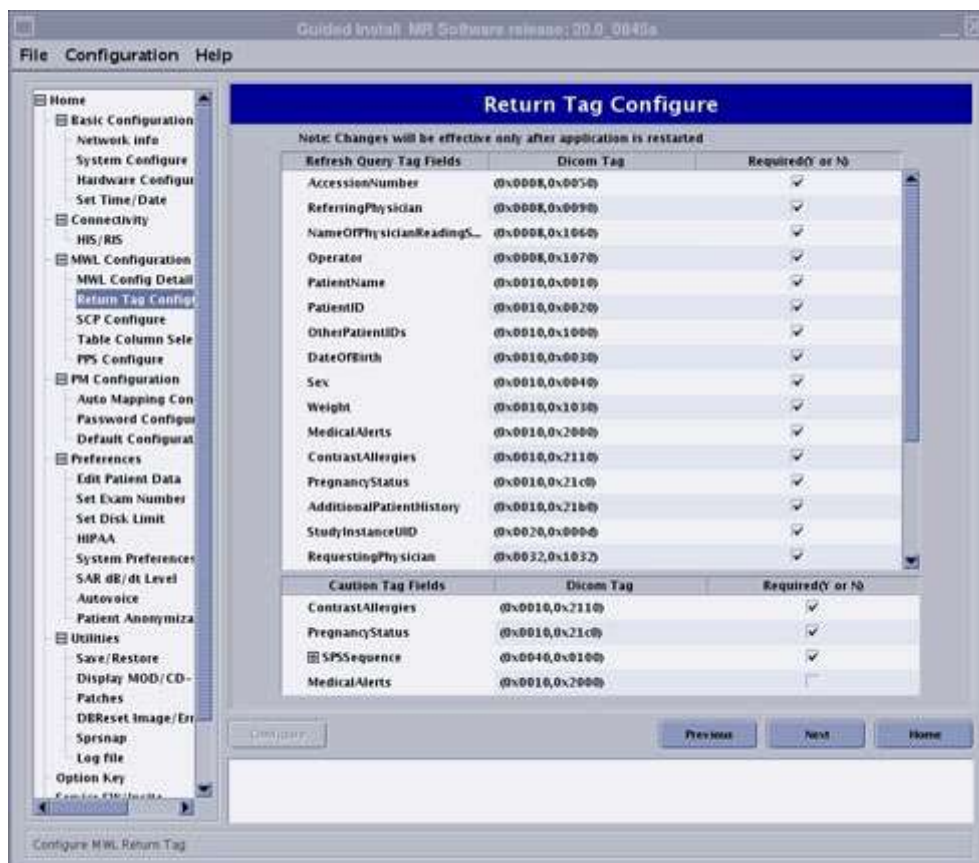

7. If any fields, other than the defaults, are desired by the site to be refreshed in the Modality Worklist, select or deselect appropriately and click the [Configure] button.
8. Under [MWL Configuration], select [SCP Configure].

Illustration 5: SCP Configure Screen

Guided Install MWL Software release: 20.0.0040a

File Configuration Help

Home

- Basic Configuration
 - Network Info
 - System Configure
 - Hardware Configure
 - Set Time/Date
 - Connectivity
 - HIS/RIS
 - MWL Configuration
 - MWL Config Detail
 - Return Tag Configure
 - SCP Configure
 - Table Column Selection
 - PPS Configure
 - PM Configuration
 - Auto Mapping Configuration
 - Password Configuration
 - Default Configuration
 - Preferences
 - Edit Patient Data
 - Set Exam Number
 - Set Disk Limit
 - HIPAA
 - System Preferences
 - SAR dR/dt Level
 - Autovoice
 - Patient Anonymization
 - Utilities
 - Save/Restore
 - Display MOD/CD-Patches
 - DEReset Image/Eraser
 - Sprsnap
 - Log file
 - Option Key

SCP Configure

Host Name	IP Address	Port Number	AE Title	DefaultSelected
RisServer	3.7.27.171	6104	RisServer	☒

Host Name: IP Address:

Port Number: AE Title:

☐ DefaultSelected


New Edit Find Remove Cancel

Previous Next Home

SCP Configure Panel

9. Click the [New] button to setup a new item.
10. Enter applicable information in all fields (as provided by site) and click the [Configure] button.
11. To verify network connectivity, click the [Ping] button.
12. Under [MWL Configuration], select [Table Column Selection Configure].

Illustration 6: Table Column Selection Configuration Screen

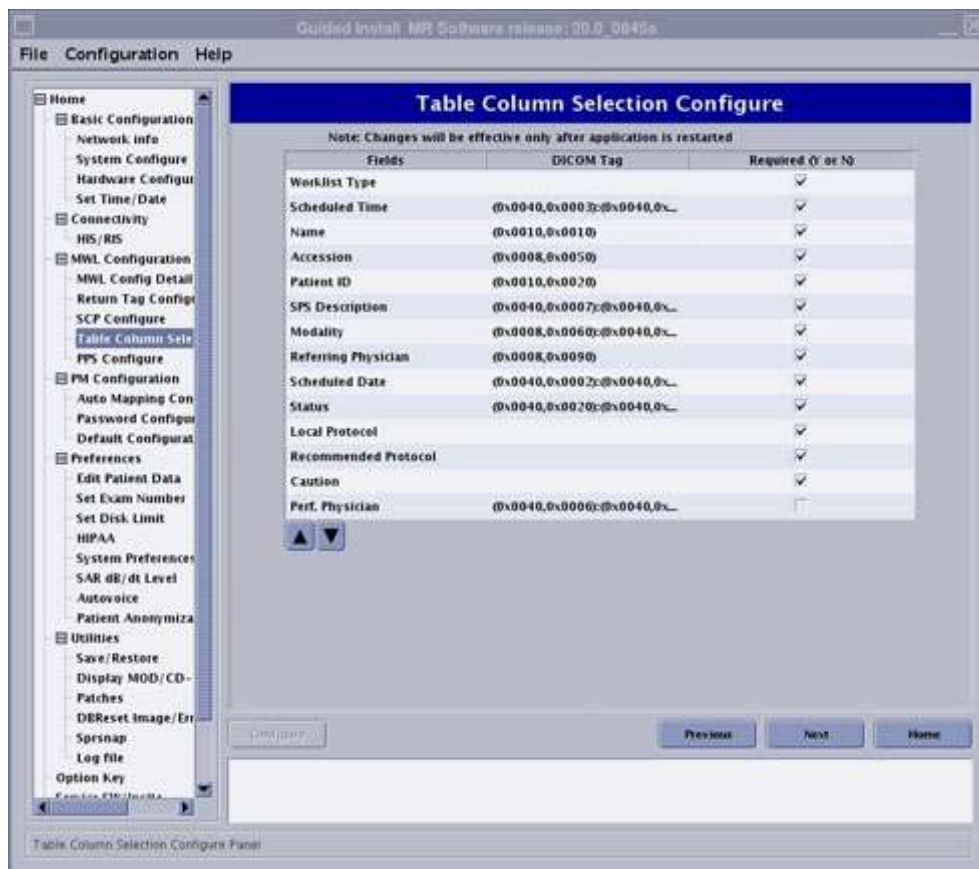


13. Select or deselect any fields desired to be displayed in the patient record table (located at the top of the scan window) and click the [Configure] button.

14. Under [MWL Configuration], select [PPS Configure].

NOTE: At release of this document, PPS Configure is available with 1.5T systems only. PPS Configure will be available for 3.0T at a later time.

Illustration 7: PPS Configure Screen

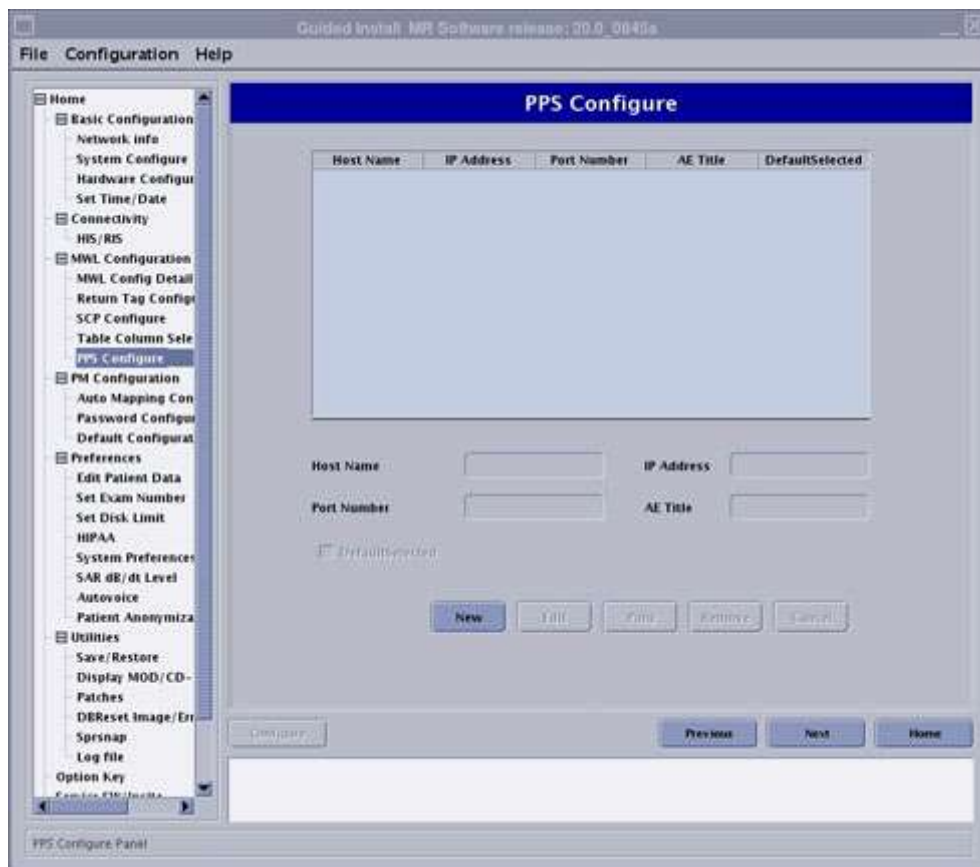


15. Click the [New] button to setup a new setup.

16. Enter applicable information in all fields (as provided by site) and click the [Configure] button.

4.2 Functional Check

1. Click the [Work List] icon for the scan screen (see [Illustration 8](#)).

Illustration 8: Work List Icon

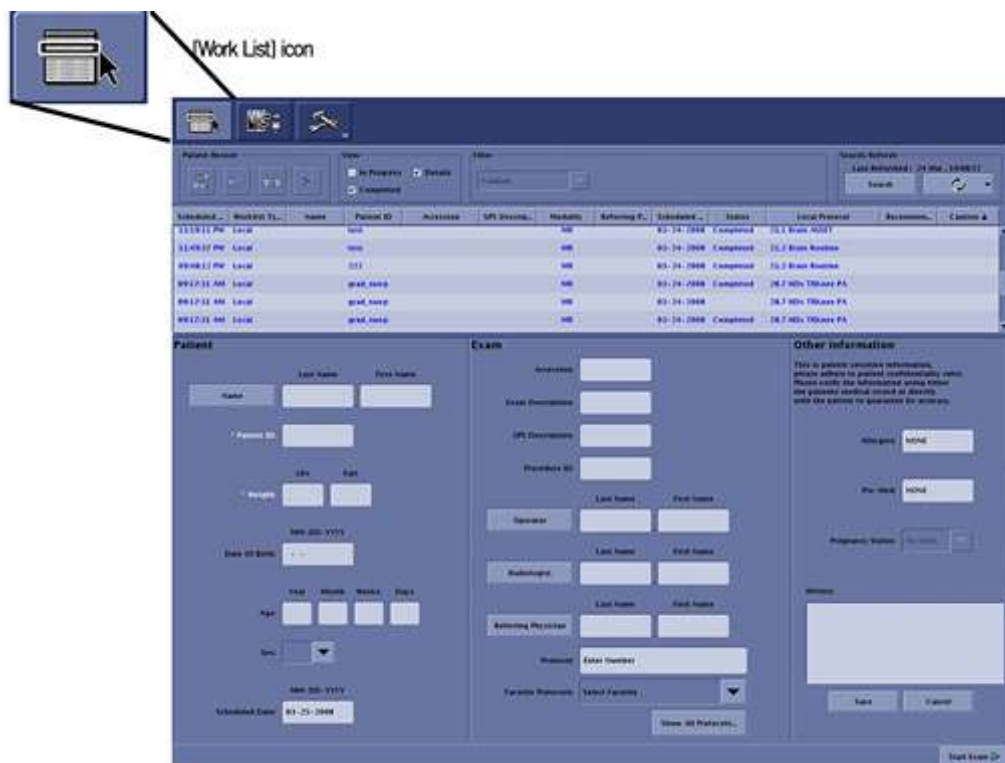


2. To confirm the connectivity of the system to the RIS system, go into the scan desktop and select the [New Patient] icon (see [Illustration 9](#)).

Illustration 9: New Patient Icon



3. On the scan screen enter the following:
 - a. **Patient ID:** *geservice*
 - b. **Weight:** *111*
4. Click the [Save] button.
5. In the Search/Refresh box, click the [Search Data] button (see [Illustration 10](#)).

Illustration 10: Search Button

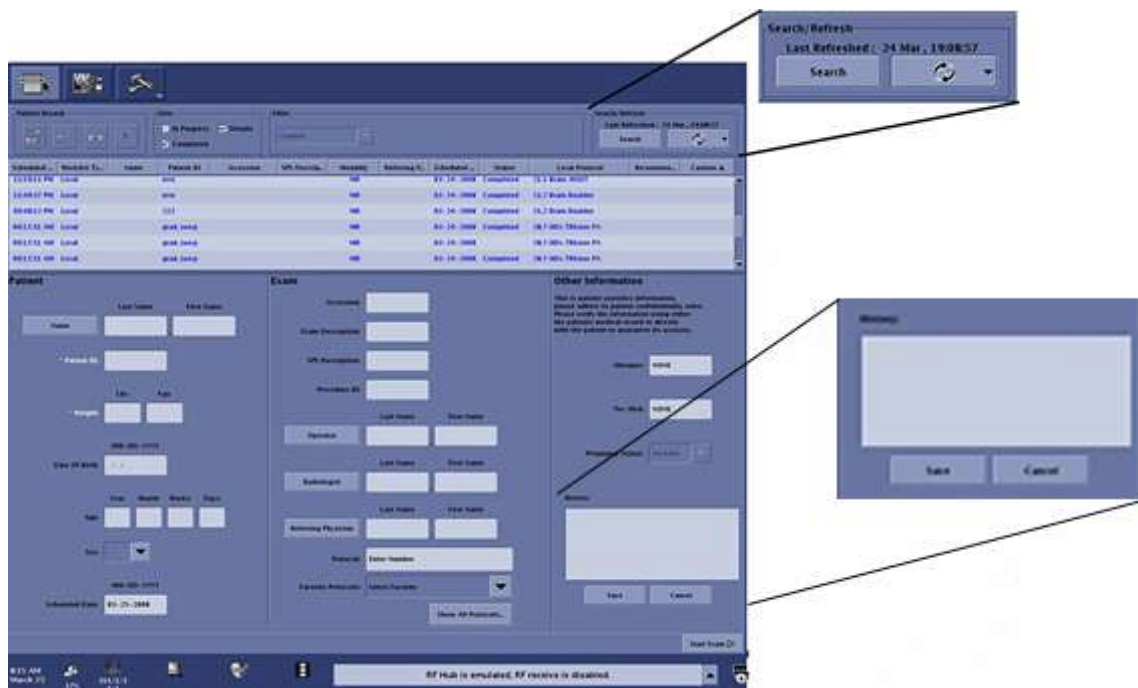

The Search Data window displays as shown in [Illustration 11](#).

Illustration 11: Search Data Window

The screenshot shows a window titled "Search Data". It contains several input fields and a section for search criteria. At the top, there are four input fields labeled "Name", "Patient ID:", "Accession:", and "Procedure ID:". Below these is a section titled "Breadth of Search". This section is divided into two columns. The left column is titled "With Date Range:" and contains five radio button options: "Today" (selected), "Any Day", "Days Before today" (with a numeric input field set to 0), "Days After Today" (with a numeric input field set to 0), and "Date Range". The "Date Range" option has "From" and "To" date input fields, both with a placeholder "MM-DD-YYYY". The right column is titled "Search RIS for:" and contains three radio button options: "This System" (selected), "This Modality" (with a text input field and the label "Enter AE Title (Optional)"), and "All Modalities". At the bottom of the window are three buttons: "Clear", "Search", and "Close".

From this page select how you want to pull information from the RIS system. The Get Patient List options are:

1. Enter a patient Name, Patient ID, Accession Number, or Procedure ID to define the search.
2. Selecting This System, which pulls the patient lists entered for this system.
3. Selecting This Modality pulls the patient list for all MR systems at this facility.
4. Selecting All Modalities pulls the patient list for all systems associated with this facility. Using this option may cause the scanner to hang if the facility has too large of Work List for all scanners at the site.
5. Selecting Date Range pulls patient lists for any set day or date range.
6. After entering the proper Search Data, click the [Search] button (see [Illustration 12](#)).

Illustration 12: Search Details Selected

Search Data

Name:

Patient ID:

Accession:

Procedure ID:

Breadth of Search

With Date Range

☒ Today

☐ Any Day

☐ Days Before today:

☐ Days After Today:

☐ Date Range

MM-DD-YYYY

From:

MM-DD-YYYY

To:

Search RIS for:

☒ This System

☐ This Modality:

Enter AE Title (Optional)

☐ All Modalities

Clear Search Close

7. The results appear in the list above the patient information (see [Illustration 13](#)).

Illustration 13: Full Work List


Search Results

Scheduled	Modality	Accession	Procedure ID	Accession	APR Study	Modality	Referring C.	Scheduled	Status	Local Patient	Accession	Accession
09:00:00 PM	Local	MRECHO	1000	1000	1000	1000	1000	09:00:00 PM	Completed	24.8 MB - 0.0s		
09:00:00 PM	Local	MRECHO	1000	1000	1000	1000	1000	09:00:00 PM	Completed	24.8 MB - 0.0s		
09:00:00 PM	Local	MRECHO	1000	1000	1000	1000	1000	09:00:00 PM	Completed	24.8 MB - 0.0s		
09:00:00 PM	Local	MRECHO	1000	1000	1000	1000	1000	09:00:00 PM	Completed	24.8 MB - 0.0s		
09:00:00 PM	Local	MRECHO	1000	1000	1000	1000	1000	09:00:00 PM	Completed	24.8 MB - 0.0s		
09:00:00 PM	Local	MRECHO	1000	1000	1000	1000	1000	09:00:00 PM	Completed	24.8 MB - 0.0s		

4.3 Troubleshooting


1. From the Common Service Desktop, select the [Guided Install] button, and in FE Mode login as *root* (see [Illustration 14](#)).

Illustration 14: Accessing HIS/RIS Setup

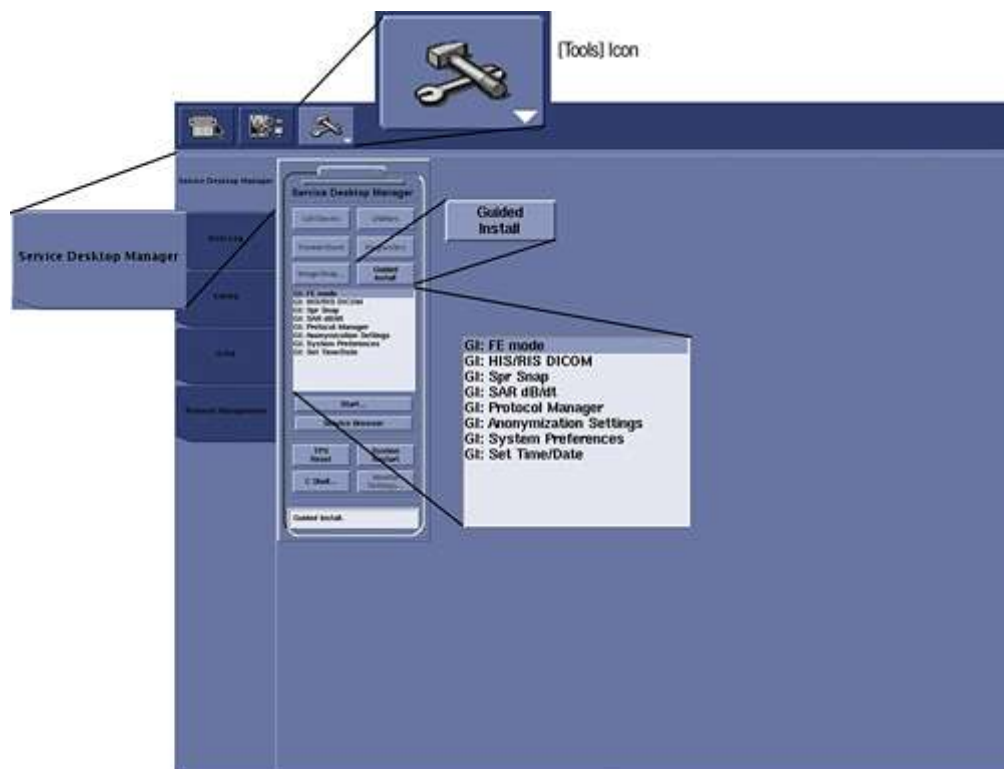

2. Under [Connectivity] select [HIS/RIS] (see [Illustration 15](#)).


Illustration 15: HIS/RIS Setup Screen and Tabs

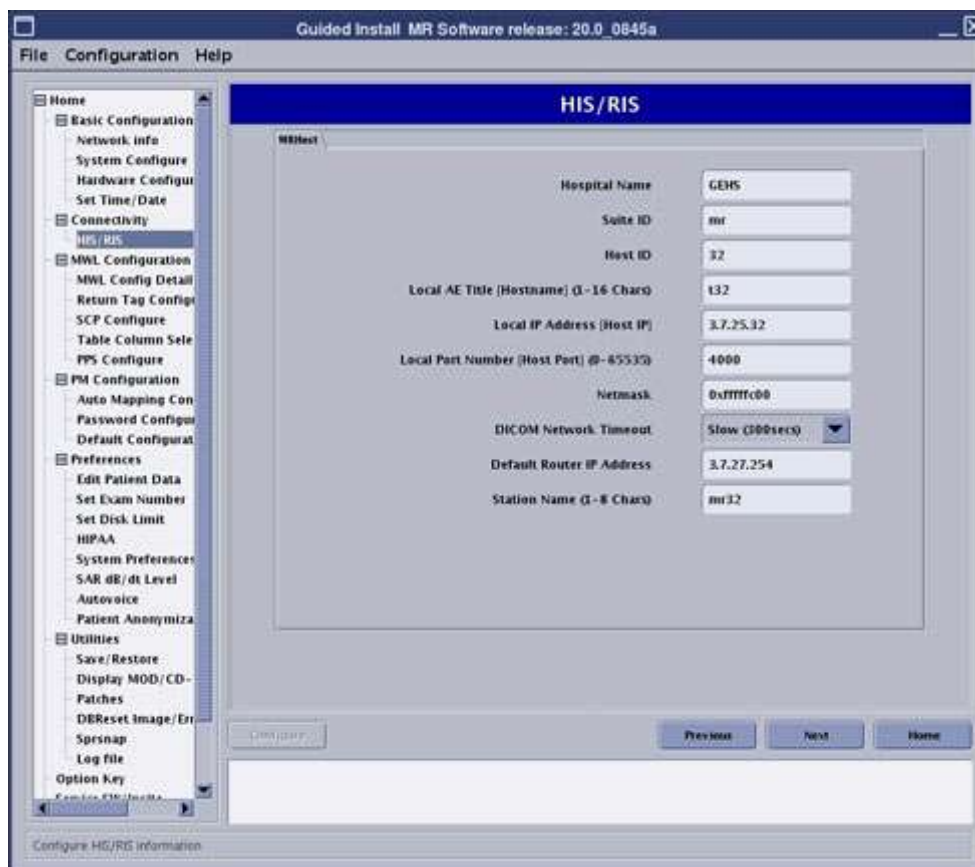

3. Ensure that the **Hospital Name**, **Suite ID**, **Host ID**, and **Local AE Title** are correct for the MR system. The **Local Port Number (Host Port)** should be set to **4000** (same as the DICOM (Digital Imaging and Communications in Medicine) networking port number for MR systems). These are the parameters that the hospital RIS system uses to identify our system (see [Illustration 15](#)).

4. Under [MWL Configuration], select [SCP Configuration] (see [Illustration 16](#)).

Illustration 16: SCP Configure Screen

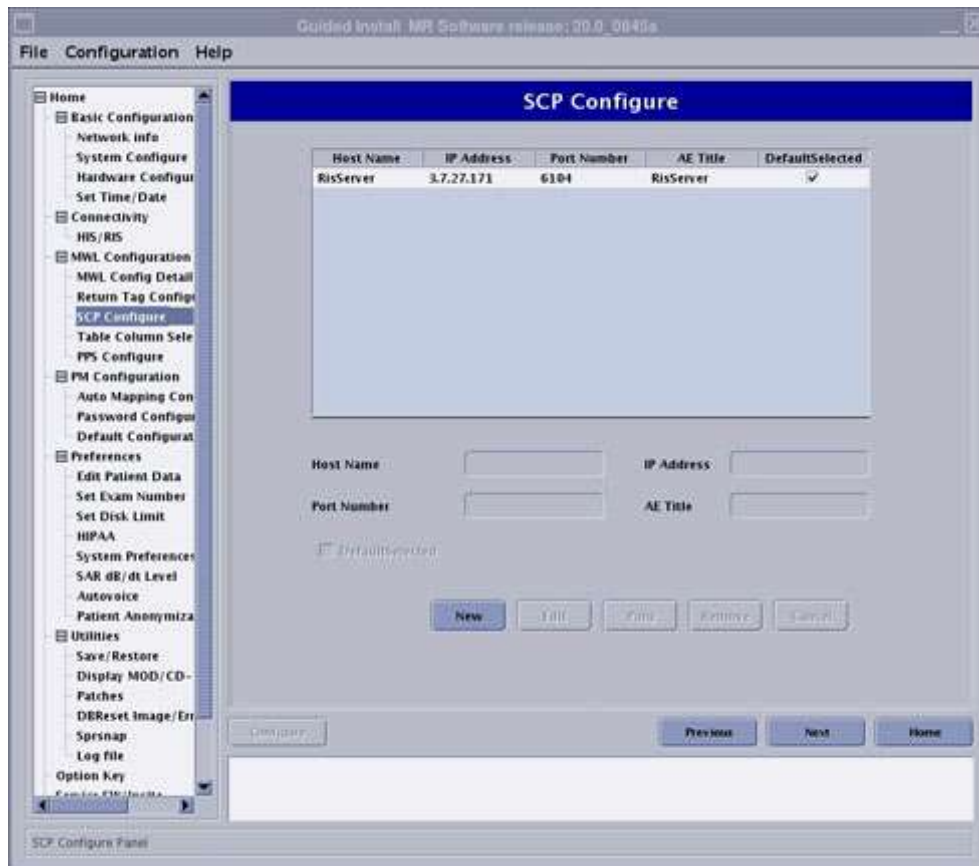


Ensure the following information with the site network administrator:

- **Host Name** is the site name.
- **AE Title** is the Modality Work List (MWL) Application Entity (AE) title of the RIS system of the site.
- **IP Address** is the internet protocol (IP) address of the RIS server.
- **Port Number** is the port number of the RIS system of the site.

5. Under [Home], select [Verification] (see [Illustration 17](#)).

Illustration 17: Verification Screen

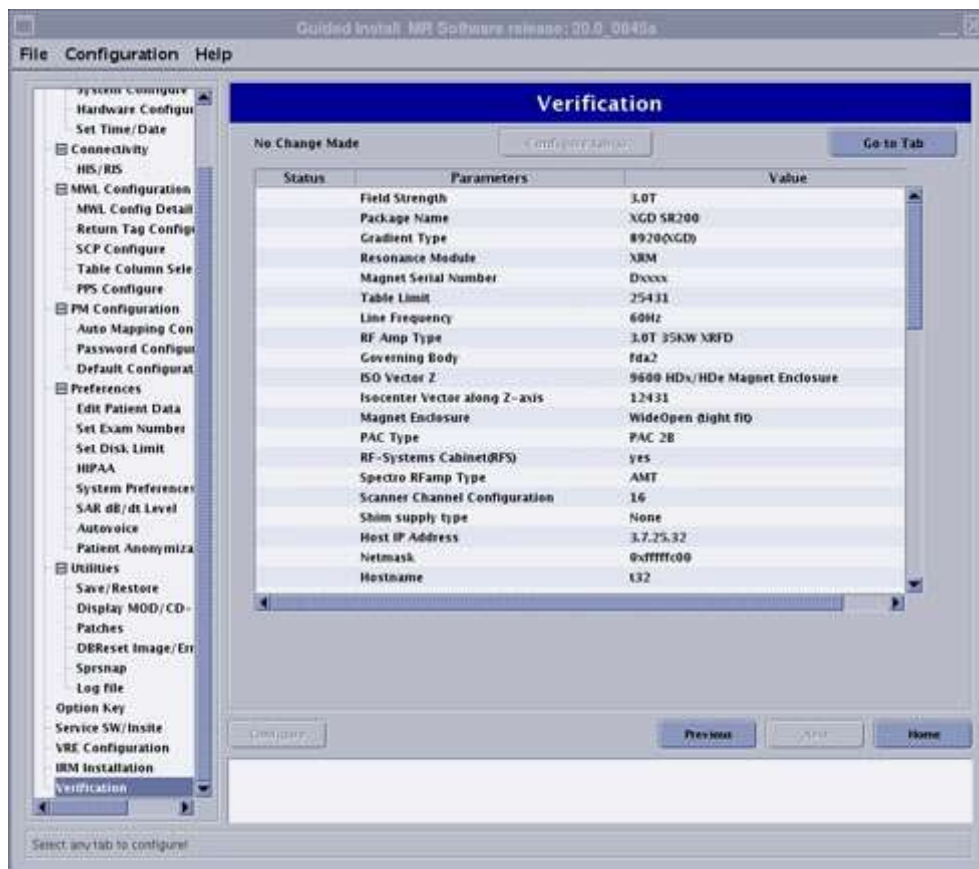

6. Select [Configure] and confirm no errors.

7. Select File>Quit.

5 Finalization

No finalization steps.

 Back to top