





1.5T SIGNA TWINSPEED – W/EXCITE HD

RECOMMENDED MAGNET ROOM
CEILING HEIGHT 9'-1" [2768mm]

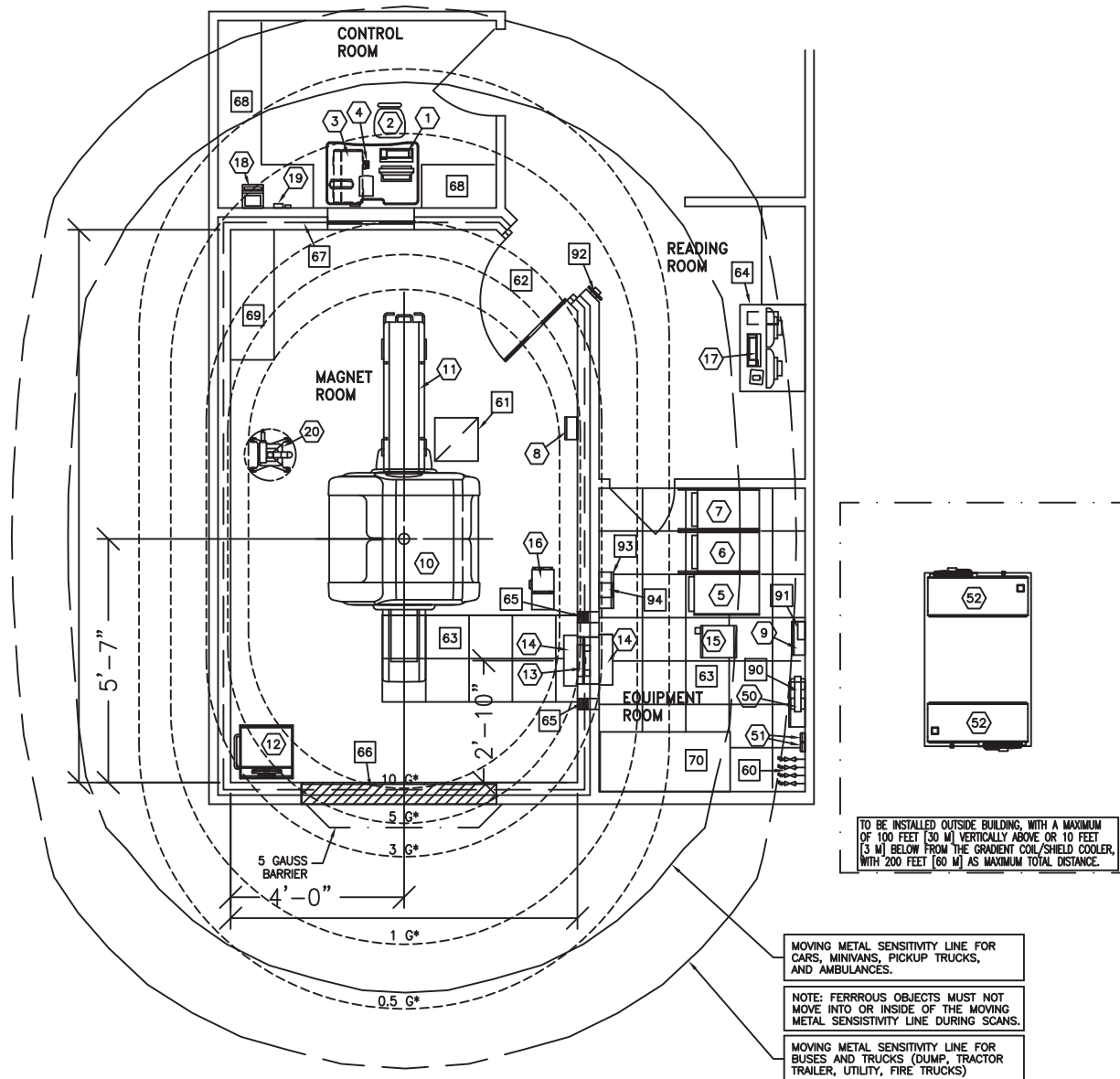

ROOM DIMENSIONS	LENGTH & WIDTH	CEILING HEIGHT
MIN. MAGNET RM:	19'-8" X 11'-0" [5.99m X 3.35m]	8'-9" [2.67m]
MIN. EQUIP. RM:	8'-0" X 12'-0" [2.44m X 3.67m]	8'-0" [2.44m]

REFER TO PRE-INSTALLATION MANUAL FOR ROOM SIZES BASED UPON CONFIGURATION

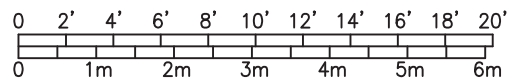

SEE APPENDIX "A" FOR POWER SPECIFICATIONS

ROOM PLAN LOCATED ON PAGE 8-189

ITEM ⬡	QUANTITY		ITEM DESCRIPTION	WEIGHT	HEAT OUTPUT
	↓				
①	1		OPERATOR WORKSPACE W/COLOR LCD MONITOR	176 lbs	4948 btu
②	1		OPERATOR'S CHAIR		
③	1		OPERATOR WORKSPACE CABINET	198 lbs	
④	1		PATIENT ALERT CONTROL BOX		
⑤	1		TWINSPEED ACCESSORY CABINET	601 lbs	2354 btu
⑥	1		RFS CABINET	780 lbs	27071 btu
⑦	1		HFD/PDU CABINET	1805 lbs	34129 btu
⑧	1		MAGNET RUNDOWN UNIT	8 lbs	
⑨	1		MAGNET MONITOR	22 lbs	204 btu
⑩	1		1.5 TESLA LCC ACTIVE SHIELD MAGNET	12826 lbs	8191 btu
⑪	1		PATIENT TRANSPORT TABLE <DOES NOT INCLUDE PATIENT>	279 lbs	
⑫	1		SPT PHANTOM CABINET	350 lbs	
⑬	1		RF PENETRATION PANEL	88 lbs	324 btu
⑭	2		PENETRATION PANEL COVERS		
⑮	1		SHIELD COOLER CABINET	275 lbs	
⑯	1		BLOWER BOX	19 lbs	1365 btu
⑰	0		ADVANTAGE WORKSTATION WITH TWO LCD MONITORS <OPTION>	81 lbs	1109 btu
⑱	0		CONTROL ROOM UNIT <OPTION>	15 lbs	
⑲	0		BATTERY CHARGING UNIT <OPTION>	4 lbs	
⑳	0		INJECTOR HEAD ON PEDESTAL <OPTION>	59 lbs	
㉔	1		MAIN DISCONNECT CONTROL	260 lbs	901 btu
㉕	2		REMOTE CONTROL FOR CHILLER SYSTEM	2 lbs	
㉖	2		MR COMMON CHILLER SYSTEM	683 lbs	

ENVIRONMENT

Ambient Operation Temperature: 50-70° F [10-21° C]
 Allowable Temperature Change: 5° F [3° C] per hour
 Humidity: 30-60 per cent non condensing
 Allowable Humidity Change: 10 per cent per hour
 Max - Altitude: 8000 ft. [2438m] above sea level

ROOM PLAN LOCATED ON PAGE 8-189

ITEM	QUANTITY	
		ITEM DESCRIPTION
		***** CUSTOMER/CONTRACTOR SUPPLIED AND INSTALLED ITEMS *****
		VALVES AND HOSE BARBS FOR COOLING SYSTEM
		MAGNET ROOM EXHAUST FAN
		MINIMUM DOOR OPENING FOR EQUIPMENT DELIVERY IS 43 IN. W x 82 IN. H [1092mm x 2083mm], CONTINGENT ON A 96 IN. [2438mm] CORRIDOR WIDTH
		NON-METAL ACCESS FLOOR WITH 2' x 2' (610 x 610mm) REMOVABLE PANELS, MIN. 0'-8" (203mm) DEEP. NON-METAL ACCESS FLOORING (PANELS & SUPPORT HARDWARE) REQUIRED WITHIN MAGNET ROOM.
		WORKSTATION TABLE
		RF FILTERS - LOCATE WITHIN 24 in. [610 mm] OF THE PENETRATION PANEL.
		MINIMUM 9 FT. -0 IN. [2743 mm] x 9 FT. -0 IN. [2743 mm] REMOVABLE WALL SECTION FOR MAGNET DELIVERY/REMOVAL.
		RF SCREEN, INCLUSIVE OF WALLS, FLOOR, DOOR, ETC. GROUND IMPEDANCE GREATER THAN 1000 OHMS. ATTENUATION 100dB AT 10-150MHZ PLANEWAVE.
		COUNTERTOP WITH DRAWERS FOR MISCELLANEOUS ITEMS.
		BASE CABINET FOR STORAGE OF: SURFACE COILS, PATIENT POSITIONING PADS, PHANTOMS, ETC.
		AIR CONDITIONING. (VIBRATION ISOLATION IS RECOMMENDED AT SUPPORTS OF EACH UNIT TO BE INSTALLED.)
		MAGNET MONITOR UPS GE CAT. NO. E4504AG 50 LBS. , 450 BTU.
		MULTIPLEXER BOX (MUX)
		METAL DETECTOR (HAND HELD)
		DC LIGHTING CONTROL PANEL 155 lbs (70 kg) 1024 BTU/HR. (CAT. NO. E4502SC/SE - BASIC SYSTEM)
		DC LIGHTING AUTO TRANSFORMER 60 lbs [27 kg] (PART OF VARIABLE DIMMER SYSTEM) (CAT. NO. E4502SD/SF INCLUDES BASIC SYSTEM)