





0.35T Signa Ovation 5

RECOMMENDED MAGNET ROOM
CEILING HEIGHT 8'-6" [2.59m]

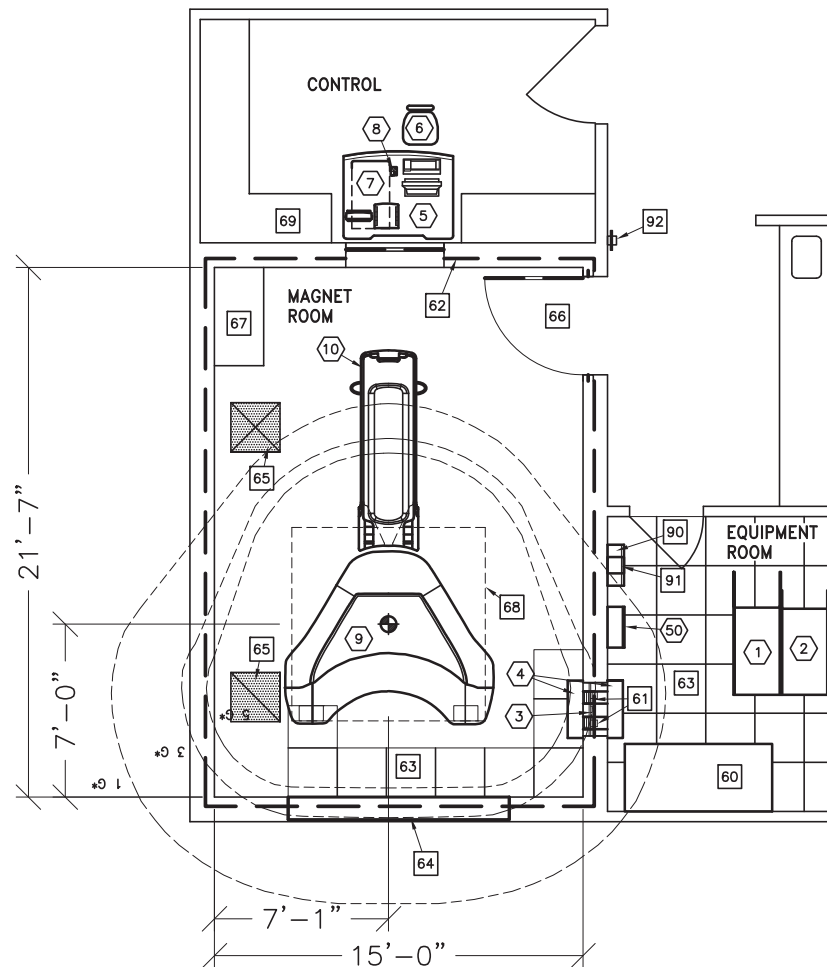

ROOM DIMENSIONS	LENGTH & WIDTH	CEILING HEIGHT
MIN. MAGNET RM:	18'-7" X 13'-1" [5.66m X 3.99m]	8'-0" [2.44m]
MIN. EQUIP. RM:	8'-0" X 7'-0" [2.44m X 2.13m]	8'-0" [2.44m]

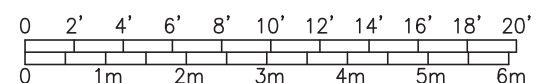

SEE APPENDIX "A" FOR POWER SPECIFICATIONS

ROOM PLAN LOCATED ON PAGE 8-217

ITEM ⬡	QUANTITY		ITEM DESCRIPTION	WEIGHT	HEAT OUTPUT
	↓				
①	1		SYSTEM CABINET	771 lbs	2542 btu
②	1		OPENSPEED POWER CABINET	1102 lbs	8866 btu
③	1		RF PENETRATION PANEL	88 lbs	324 btu
④	2		PENETRATION PANEL COVERS		
⑤	1		OPERATOR WORKSPACE W/COLOR LCD MONITOR	176 lbs	1481 btu
⑥	1		OPERATOR'S CHAIR		
⑦	1		OPERATOR WORKSPACE CABINET	191 lbs	
⑧	1		PATIENT ALERT CONTROL BOX		
⑨	1		0.35 TESLA MAGNET	41887 lbs	4709 btu
⑩	1		PATIENT TABLE	1302 lbs	

ENVIRONMENT

Ambient Operation Temperature: 59-89° F [15-32° C] – EQUIPMENT ROOM
 Allowable Temperature Change: 5° F [3° C] per hour
 Humidity: 30-60 per cent non condensing
 Allowable Humidity Change: 10 per cent per hour
 Max – Altitude: 11,808 ft. [3600m] above sea level

ROOM PLAN LOCATED ON PAGE 8-217-177

ITEM 	QUANTITY 	ITEM DESCRIPTION	WEIGHT	HEAT OUTPUT
		***** CUSTOMER/CONTRACTOR SUPPLIED AND INSTALLED ITEMS ***** [60] AIR CONDITIONING UNIT BY OTHERS. (VIBRATION ISOLATION IS RECOMMENDED AT FLOOR CONNECTION POINTS OF AIR CONDITIONING UNITS TO BE INSTALLED INSTALLED FOR COOLING THE MR SUITE. > [61] RF FILTERS BY OTHERS - LOCATE WITHIN 2'-0" [610mm] OF THE PENETRATION PANEL. [62] RF SCREEN, INCLUSIVE OF WALLS, DOOR, FLOOR, ETC. GROUND IMPEDANCE GREATER THAN 1000 OHMS. ATTENUATION 100DB AT 10-100MHZ PLANEWAVE. [63] ACCESS FLOORING BY OTHERS [64] MINIMUM 9 FT. [2743 mm] x 8 FT. [2438 mm] REMOVABLE WALL SECTION. FOR MAGNET DELIVERY/REMOVAL. [65] AIR SUPPLY AND RETURN DUCTS. [66] MINIMUM DOOR OPENING FOR EQUIPMENT DELIVERY IS 43 IN. W X 82 IN. H [1092MM X 2083MM], CONTINGENT ON A 96 IN. [2438MM] CORRIDOR WIDTH [67] BASE CABINET FOR STORAGE OF: PATIENT POSITIONING PADS, SURFACE COILS, ETC. [68] SHIM PLATE [69] COUNTERTOP WITH DRAWERS [90] DC LIGHTING AUTO TRANSFORMER 60 lbs [27 kg] (USED ONLY IN VARIABLE DIMMER SYSTEM) (CAT. NO. R4503AF/AY INCLUDES CONTROL PANEL) [91] DC LIGHTING CONTROL PANEL - 155 lbs [70 kg] 1024 BTU/HR. (CAT. NO. R4503AD/AW)		