

MAGNETOM

# Startup System

Start-up / Magnet System

Aera  
Skyra

## Document Version

Siemens reserves the right to change its products and services at any time.

In addition, manuals are subject to change without notice. The hardcopy documents correspond to the version at the time of system delivery and/or printout. Versions to hardcopy documentation are not automatically distributed.

Please contact your local Siemens office to order a current version or refer to our website <http://www.healthcare.siemens.com>.

## Disclaimer

Siemens provides this documentation “as is” without the assumption of any liability under any theory of law.

The content described herein requires superior understanding of our equipment and may only be performed by qualified personnel who are specially trained for such installation and/or service.

## Copyright

“© Siemens, 2010” refers to the copyright of a Siemens entity such as:

Siemens Healthcare GmbH - Germany

Siemens Aktiengesellschaft - Germany

Siemens Shenzhen Magnetic Resonance Ltd. - China

Siemens Shanghai Medical Equipment Ltd. - China

Siemens Healthcare Private Ltd. - India

Siemens Medical Solutions USA Inc. - USA

Siemens Healthcare Diagnostics Inc. - USA and/or

Siemens Healthcare Diagnostics Products GmbH - Germany

|          |                                                            |           |
|----------|------------------------------------------------------------|-----------|
| <b>1</b> | <b>Verifying the Electronic Low Pressure Switch</b>        | <b>5</b>  |
| <b>2</b> | <b>Filling the Magnet</b>                                  | <b>6</b>  |
| 2.1      | Quench tube check                                          | 6         |
| 2.2      | Filling the static system                                  | 7         |
| 2.2.1    | General                                                    | 7         |
| 2.2.2    | Safety                                                     | 9         |
| 2.2.3    | Site preparation                                           | 12        |
| 2.2.4    | System preparation                                         | 12        |
| 2.2.5    | Assemble the service syphon                                | 20        |
| 2.2.6    | Fit the magnet service syphon                              | 21        |
| 2.2.7    | Precool the service syphon                                 | 25        |
| 2.2.8    | Liquid helium transfer - top fill                          | 27        |
| 2.2.9    | Removing the service syphon                                | 33        |
| 2.2.10   | Changing supply dewars during transfer                     | 36        |
| 2.2.11   | Stopping the helium transfer                               | 39        |
| 2.2.12   | Final work steps                                           | 40        |
| 2.2.13   | Low ceiling height                                         | 41        |
| 2.2.14   | Helium pressures at altitude                               | 42        |
| 2.2.15   | Verifying the electronic low pressure switch               | 44        |
| 2.3      | Bottom filling the magnet                                  | 50        |
| 2.3.1    | Checking magnet temperatures                               | 51        |
| 2.3.2    | Bottom fill procedure                                      | 54        |
| 2.3.3    | Liquid helium transfer - bottom fill                       | 57        |
| 2.3.4    | Stopping the helium transfer                               | 63        |
| 2.3.5    | Removing the service syphon                                | 63        |
| 2.3.6    | Final work steps                                           | 66        |
| <b>3</b> | <b>Prerequisites for Ramping the Magnet</b>                | <b>68</b> |
| 3.1      | Quench tube check                                          | 68        |
| 3.2      | Identification of the controlled access area (0.5 mT area) | 69        |
| 3.3      | Magnet stop button test                                    | 70        |
| 3.4      | Venting leak check                                         | 71        |
| 3.4.1    | General                                                    | 71        |
| 3.4.2    | Venting System                                             | 71        |
| 3.4.3    | Safety                                                     | 72        |
| 3.4.4    | Tooling and Spare Parts                                    | 73        |
| 3.4.5    | Leak Check Process at Installation and Service             | 74        |
| 3.4.6    | High Pressure Turret and Venting Components                | 76        |
| 3.4.7    | Internal Magnet Leak Checks                                | 84        |
| 3.4.8    | Low Pressure Turret and Venting Components                 | 90        |
| 3.4.9    | Final Work Steps                                           | 92        |
| <b>4</b> | <b>Ramping the Magnet</b>                                  | <b>93</b> |
| 4.1      | OR98/99/103 Ramping                                        | 93        |
| 4.1.1    | Tools required                                             | 95        |
| 4.1.2    | Preparation                                                | 96        |
| 4.1.3    | Installation                                               | 97        |

|          |                                                          |            |
|----------|----------------------------------------------------------|------------|
| 4.1.4    | Set up of the MPSU .....                                 | 100        |
| 4.1.5    | Run to Field .....                                       | 122        |
| 4.1.6    | Run to Zero .....                                        | 126        |
| 4.1.7    | Ramping the magnet to the required frequency. ....       | 128        |
| 4.1.8    | Work after ramping the magnet .....                      | 130        |
| 4.1.9    | Final work steps. ....                                   | 131        |
| 4.1.10   | Mains failure during ramping .....                       | 134        |
| 4.1.11   | iButton function .....                                   | 135        |
| 4.2      | What to do after a quench. ....                          | 136        |
| 4.2.1    | General .....                                            | 136        |
| 4.2.2    | Post Quench Steps .....                                  | 137        |
| 4.2.3    | Pre-requisites. ....                                     | 140        |
| 4.2.4    | Checking the Quench Valve and Burst Disk. ....           | 150        |
| 4.2.5    | Checking the Auxiliary Burst Disk .....                  | 154        |
| 4.2.6    | Final Work Steps .....                                   | 158        |
| <b>5</b> | <b>Shimming the Magnet .....</b>                         | <b>163</b> |
| 5.1      | Prerequisites for Shimming .....                         | 163        |
| 5.1.1    | Only for Systems with a Tim Dockable Patient Table ..... | 163        |
| 5.2      | Passive Shim .....                                       | 165        |
| <b>6</b> | <b>Changes to Previous Version .....</b>                 | <b>166</b> |
| 6.1      | Changes to previous version .....                        | 166        |
| <b>7</b> | <b>List of Hazard IDs .....</b>                          | <b>168</b> |

Verifying the electronic low pressure switch is part of the filling procedure.

If the magnet does not require a helium fill, the test must be performed as a stand-alone procedure.



If the magnet does not require a helium fill at start-up, the Low Pressure Switch must be tested.

Refer to (→ Verifying the electronic low pressure switch / Page 44) in the Filling Procedure

## 2.1 Quench tube check

- Perform a visual inspection of the quench tube.
  - » Refer to (→ Quench tube / M7-010.833)
- Check if the Installation protocol is signed by the project manager and Quench line manufacturer.



Ensure the Quench Line installation is complete before performing the first filling procedure.

The Quench Line Acceptance Protocol must be checked before the first ramp.

Check that the Quench Line Acceptance Protocol is signed off by the Quench Line Manufacturer.

The acceptance protocol can be found under:

**Product Information Home > Planning > Planning Guides > MR Systems > Quench Line Design and Calculation Tool**

## 2.2 Filling the static system


### 2.2.1 General

#### 2.2.1.1 Static systems

The magnet design does not include a permanently fitted syphon. A short magnet service syphon must be fitted when filling with liquid helium. The service syphon is tooling only and must be removed immediately after filling has been completed.

This procedure can be performed on magnets at field and magnets at zero field.

Fig. 1: Upgrade kit 10129438 - De-ice tooling

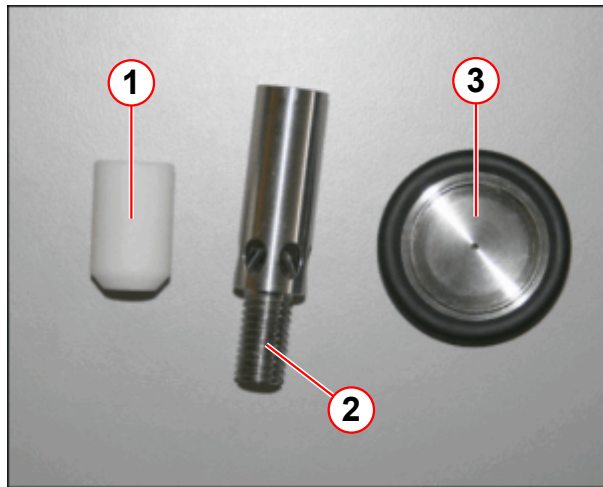

- (1) PTFE tip
- (2) De-ice helium fill diffuser
- (3) NW25 flow restrictor

The Magnet service syphon kit is now upgraded to include **de-ice tooling**.

The de-ice tooling is used for **helium top fill**, and **de-icing of the syphon port**.

- Use for **helium top filling**:
  - PTFE tip
  - 6 hole metal de-ice helium fill diffuser
- Use for **de-icing of the syphon port**:
  - PTFE tip
  - 6 hole metal de-ice helium fill diffuser
  - NW25 flow resistor

The new 6 hole diffuser is **more efficient**, and must be used **in place of** the 60 mm PTFE diffuser.



Only Siemens trained and accredited CSEs are allowed to perform this work.



During filling, do not exceed 95%.

| System Type           | Recommended helium level for the first ramp to field | Minimum helium level for future ramping | Minimum helium level for customer hand over |
|-----------------------|------------------------------------------------------|-----------------------------------------|---------------------------------------------|
| Aera OR98             | 60%                                                  | 40%                                     | 50%                                         |
| Skyra OR99            | 65%                                                  | 60%                                     | 50%                                         |
| Skyra / Spectra OR103 | 65%                                                  | 40%                                     | 50%                                         |



To recover a potentially warm magnet, refer to magnet temperatures in the Bottom filling procedure (→ Bottom filling the magnet / Page 50). Contact Magnet HSC for advice.



It is recommended that fills are performed during service visits and before ramping up.  
The minimum level for filling at field is 40%  
Wait at least 2 hours after filling before ramping the magnet.



If a helium fill is showing poor fill efficiency, or no increase on the helium level meter, the transfer must be stopped and the syphon port inspected for ICE.

#### 2.2.1.2 Venting system



#### CAUTION

Cold head failure.


The cold head seals can fail, and leak to atmosphere.

» Always ensure the cold head bypass valve is CLOSED during helium filling.



This document describes the filling procedure for **OR98, OR99 and OR103** magnets.

Fig. 2: Normal venting configuration

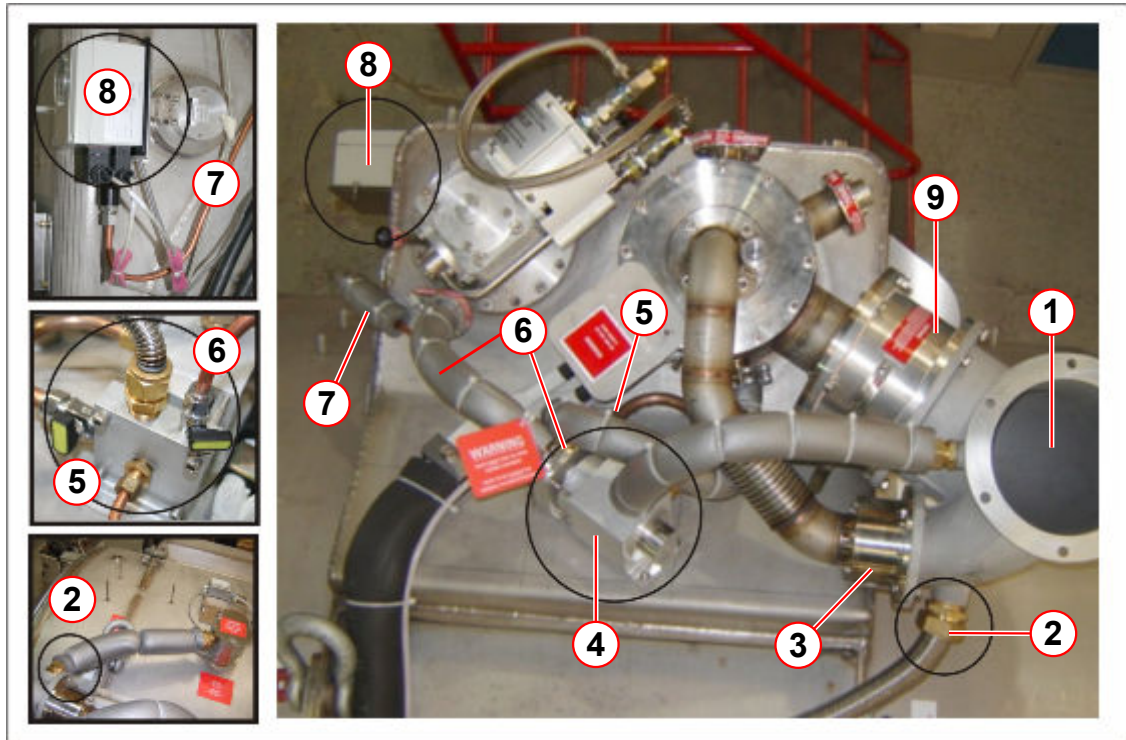

- (1) Quench valve elbow
- (2) Flexible connection from OVC vent port
- (3) Auxiliary burst disk (1.38 bar)
- (4) Venting valve assembly
- (5) Magnet system bypass line
- (6) Cold head vent line and bypass valve
- (7) Pressure display & control switch line
- (8) Pressure display & control switch
- (9) Quench valve with integral 15 psi (G) burst disk

## 2.2.2 Safety



### WARNING

**Strong magnetic field.**

**Failure to observe the following safety warning may result in physical injury and equipment damage.**

- » Only non-magnetic equipment may be taken into the magnet room.
- » Never transport the GHe steel cylinder into the magnet room.



Ensure that the Quench Line has been correctly installed and safe before performing Installation and Startup procedures.

**WARNING**

Risk of asphyxiation.

Failure to observe the following safety warning may result in dizziness and loss of consciousness in personnel working with cryogen.

- » Do not vent helium gas into the magnet room during helium transfer. As an exception, helium gas may be vented to cool the syphon.
- » Ensure that the turret venting is connected to quench pipe.
- » Oxygen monitors must be activated at all times.

**CAUTION**

Cryogenic liquids.

Failure to observe the following safety warning may result in cold burns.

- » Protective clothing, gloves and eye protection must be worn.
- » Filling must be performed by trained personnel only.

**CAUTION**

Pressurized helium vessel.

Cold helium gas.

Failure to observe the following safety information may result in physical injury.

- » Do not disconnect any of the turret venting or fittings without first depressurizing the helium vessel.

**CAUTION**

Risk of quench during filling

If the magnet quenches during helium transfer

- » Stop the helium transfer if possible and safe to do so.
- » Follow Post Quench procedures


**NOTICE**

Mishandling of the dewar.

- » The dewar-relevant safety regulations of the manufacturer must be observed.

### 2.2.2.1 Equipment required

Fig. 3: Magnet service syphon kit

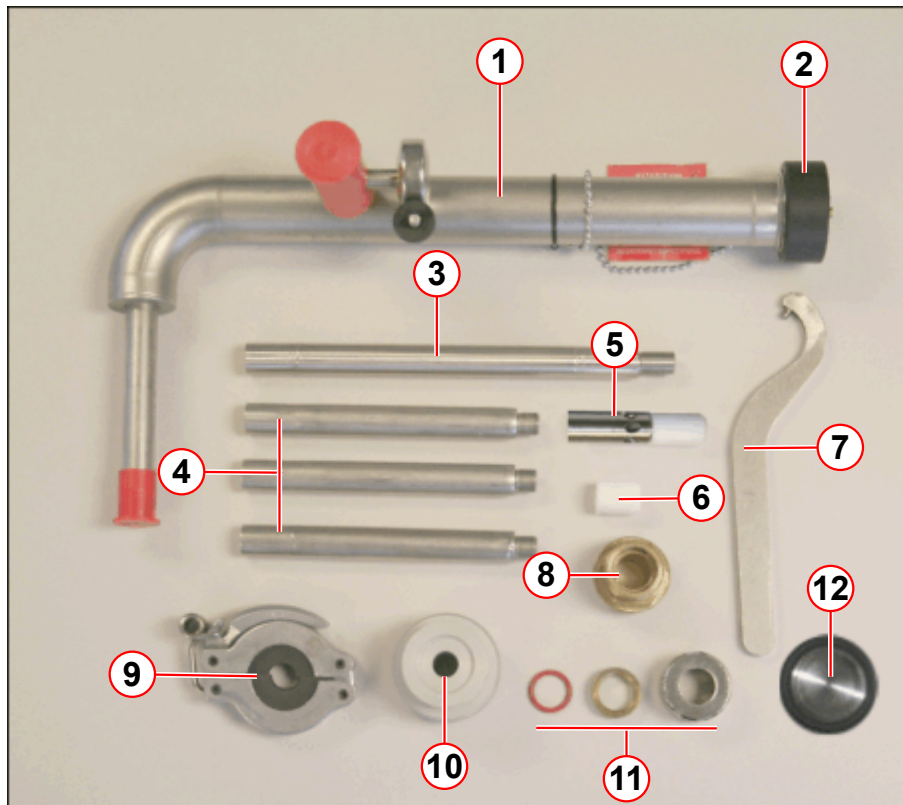

- (1) Service syphon assembly
- (2) Syphon end cap
- (3) 200 mm extension leg (1 off)
- (4) 140 mm extension legs (3 off)
- (5) 6 hole helium fill diffuser tip (top fill)
- (6) 22mm PTFE seal (bottom fill)
- (7) Syphon C-spanner/wrench
- (8) M28 brass syphon adaptor bush (OR105, OR122)
- (9) Syphon leg support clamp (low ceiling)
- (10) Nitrogen fill adaptor and O-ring seal
- (11) Syphon O-ring seal, brass syphon flange, syphon ring nut
- (12) NW25 flow restrictor

- Magnet service syphon kit (part no **10126651 (Revision 01)**). (→ Fig. 3 Page 11)
- Edwards hand held leak checker, (part number **10614850**)
- Dewar syphon.
- Dewar neck adaptor, clamps, seals, and hose.
- Heat gun.
- Paper towelling.
- Helium gas supply (at least 99.9% pure).
- Pressure regulator 200 bar/ 0 - 1.5 bar.
- Protective clothing: gloves, safety glasses, face shield, lab coat, safety shoes.
- Replacement copper seal for syphon port on the turret.

- Service platform.
- Gas bottle wrench.
- A room or personal O<sub>2</sub> Monitor.

### 2.2.3 Site preparation

- Ensure warning signs are correctly positioned.
- Ensure there is a clear path to exit the imaging suite.


### 2.2.4 System preparation

Certain operations must be performed without delay to ensure efficient transfer:-

- Ensure that there is **open and clear access** to the magnet for **helium dewar(s)**.
- Ensure that there is **open and clear access** to the **syphon**.
- **Position** the **Service platform** next to the magnet. .
  - The platform is secured to the magnet, using **2 long straps** which are passed **through the magnet bore**.
- During installation, if the **MRC (MRAWP) Host** is not available, take a helium sample using the **MSUP**. (→ Helium level check at the MSUP / Page 31)
- **Remove** the appropriate **covers**.

### 2.2.4.1 Service platform installation

Fig. 4: Service platform - as delivered

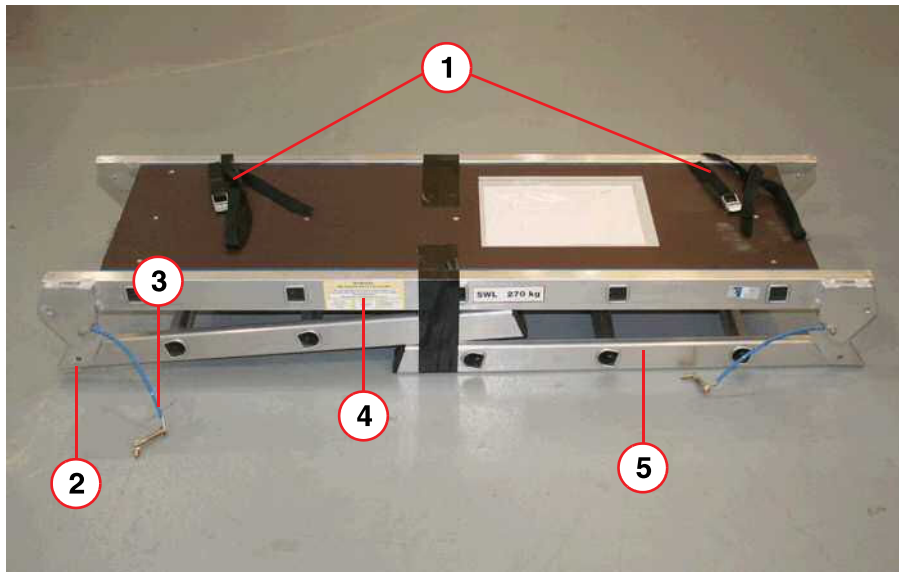

- (1) Safety straps (2 off)
- (2) Hinged locking pin location (each end)
- (3) Retaining wire and hinged locking pin (each end)
- (4) Warning label
- (5) Step legs (each end)

### Service platform safety

Before using the service platform **read the following safety cautions** and perform the **required safety checks**.



#### CAUTION

##### Risk of personal injury

Failure to observe the following precautions may result in personal injury.

- » Always check the platform and hinge locking pins for damage.
- » Always check the hinge locking pins are secured.
- » Always secure the platform to the magnet frame or with the two long straps through the magnet bore.
- » Always ascend and descend (Walk UP or DOWN) the steps facing the platform.
- » Where possible use the frame or covers for additional hand support.
- » Never climb onto the platform from the side.
- » Never jump from the steps or top of the platform.

Fig. 5: Correct method - Up to platform



Fig. 6: Incorrect method - Down from platform

**⚠ CAUTION****Risk of personal injury**

Carrying heavy components or tooling up or down the platform steps can lead to injury.


- » DO NOT carry heavy components up or down the platform steps.
- » Place heavy objects (i.e. cold head) onto the platform first.
- » Only lift the component when standing safely on the platform.

**Platform assembly**

- Remove any packaging.
- Unfold both step legs.




Fig. 7: Service platform - lock legs in position

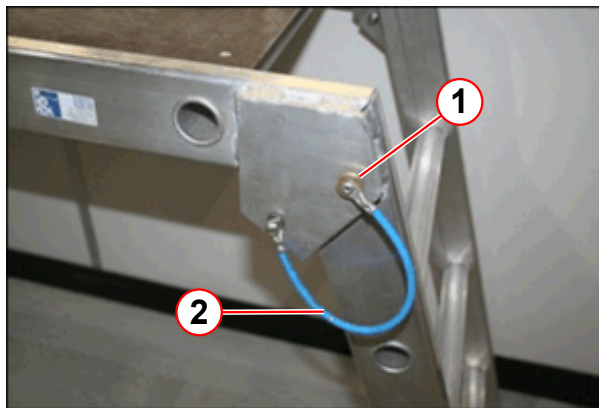

- (1) Hinged locking pin
- (2) Retaining wire


- Align the hole at the top of the leg with the **pin location hole** (pos 1).
- Push the **hinged locking pin** through the hole (pos 1) and fold back the hinge to secure.
- Repeat with the other leg.

Fig. 8: Securing the service platform (example)

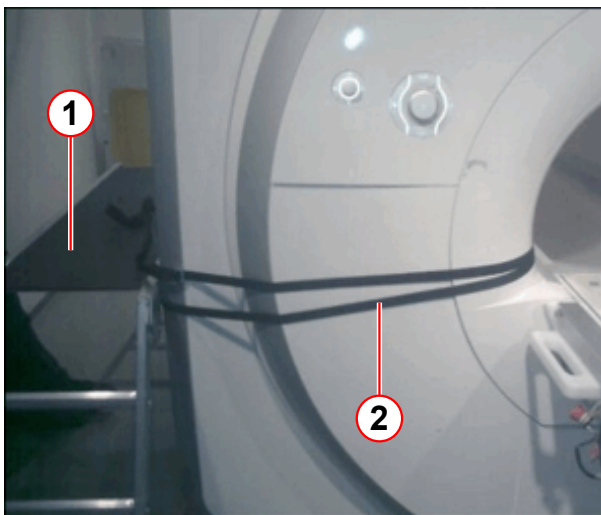

- (1) Service platform
- (2) Long straps (fitted through the magnet bore)

Fig. 9: Securing the service platform (example)

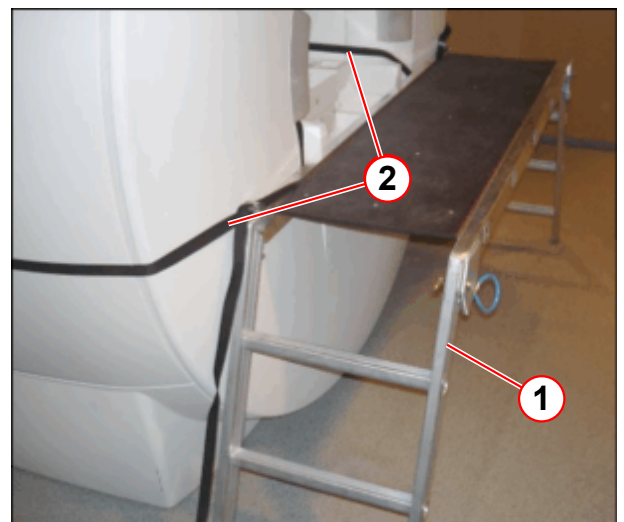

- (1) Service platform
- (2) Long straps (fitted through the magnet bore)



Do not exceed the maximum load - the combined weight of 2 people must not exceed 270 kg.



The fiber optic connections in the top of RFIS are very fragile, so exercise due caution.

#### 2.2.4.2 Dewar preparation

This procedure describes the use of **top fill** helium dewars used in conjunction with the Siemens Magnet Technology dewar syphon.

The procedure covers all dewars that are externally pressurized with helium gas.



Refer to the dewar manufacturer's documentation when using side fill and/or self-pressurizing dewars.



Refer to the dewar manufacturer's safety guidelines and check the dewar valve functions before commencing the LHe fill.


- If the dewar is **highly pressurized**, (greater than 0.3bar (4.4psi)), release the pressure by slowly opening the **primary ball valve** (→ 5/ Fig. 10 Page 16) to minimize the flash losses of helium.
  - This must be done in a **well ventilated area - NOT in the imaging suite.**



Check for the cause of the high pressure, i.e. was the low pressure relief valve closed? Inform the supplier if the cause cannot be found.

- Fit the **dewar syphon adaptor** (→ 3/ Fig. 12 Page 18) to the top of the dewar.

Fig. 10: Dewar Ports (typical)

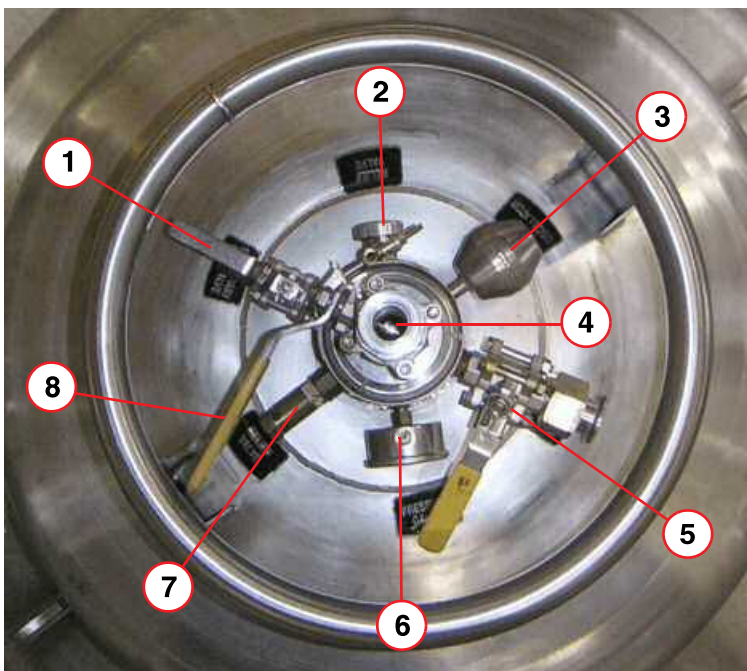


- (1) Low pressure relief valve (shown open)
- (2) Relief valve
- (3) Anti-oscillator
- (4) Syphon entry port (shown closed)
- (5) Primary ball valve (shown closed)
- (6) Pressure gauge
- (7) Relief valve
- (8) Syphon entry valve

- **Close the primary ball valve** (→ 5/ Fig. 10 Page 16) and **low pressure relief valve.** (→ 1/ Fig. 10 Page 16)



### Dewar syphon preparation

Fig. 11: Dewar syphon

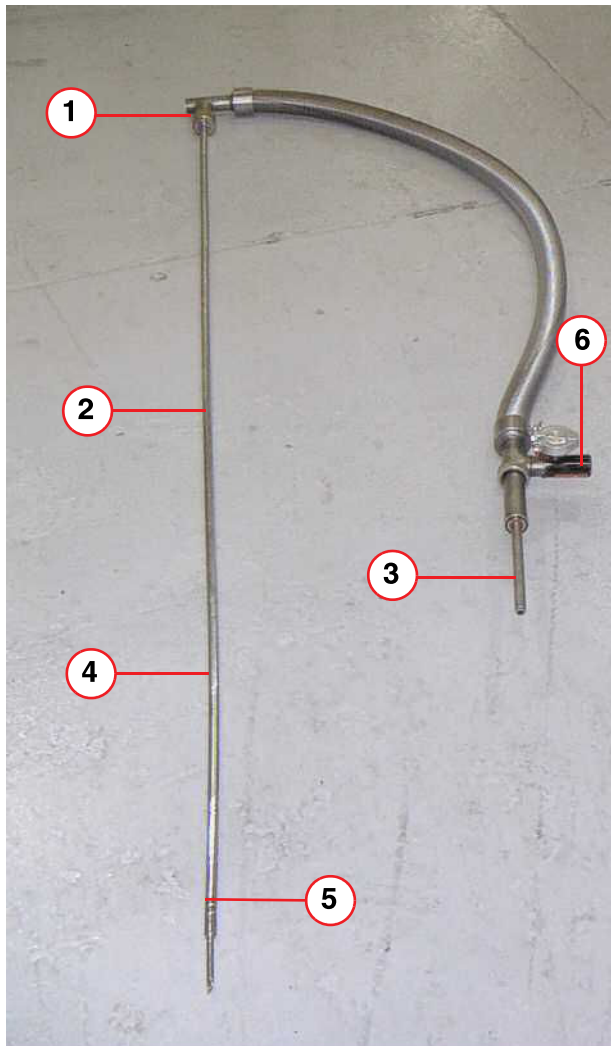

- Fit the extension leg (s) and sludge inhibitor, to the dewar **syphon leg**.

- (1) Syphon flow valve
- (2) Syphon leg
- (3) Bayonet connection (male)
- (4) Extension leg
- (5) Sludge inhibitor
- (6) Syphon vacuum valve




Do not open the syphon vacuum valve (6).

- Fit the **helium regulator** to the helium **gas bottle**.

### Precool the dewar syphon

Fig. 12: Dewar syphon fitment

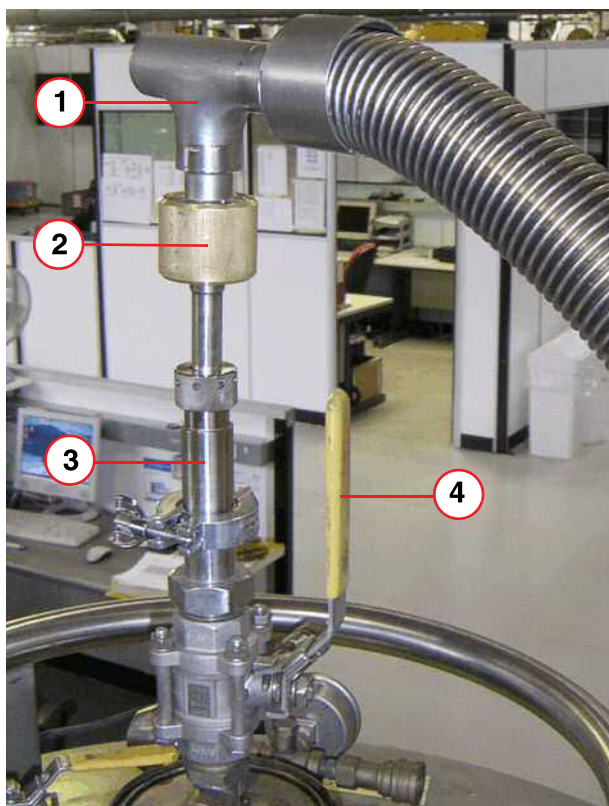


- (1) Syphon
- (2) Syphon flow valve
- (3) Dewar syphon adaptor
- (4) Syphon entry valve

- Open the **syphon flow valve** (pos 2) two turns.
- Insert the syphon leg into the **dewar syphon adaptor** (pos 3), and tighten the ring nut and seal in order to create a gas seal that the syphon can slide through.  
(Do not fully tighten the ring nut at this stage).
- Open the **syphon entry valve** (pos 4), and slowly slide the syphon into the dewar. Gas will begin to flow through the syphon.

- Push the syphon **fully down** into the dewar, (very slowly to minimize pressure build-up and flash losses), until the sludge inhibitor touches the bottom.

- Fully tighten the **ring nut and seal**.

Fig. 13: Pre-cooling the syphons

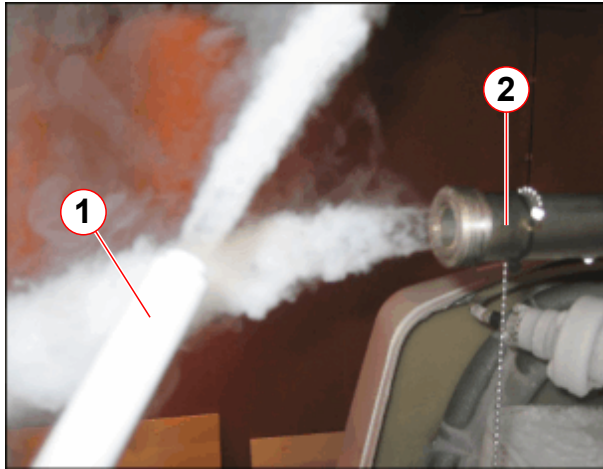


(1) Dewar syphon  
(2) Service syphon

- Allow the syphon to cool until a dense white "cascading plume" appears.
- **Do not** allow high pressure to build in the dewar, only enough to cool the syphon.

- When the syphon is cold, close the **syphon flow valve** (→ 2/ Fig. 12 Page 18) and move the dewar so that the syphon can be connected without excessive bending.

Fig. 14: Dewar gas connection

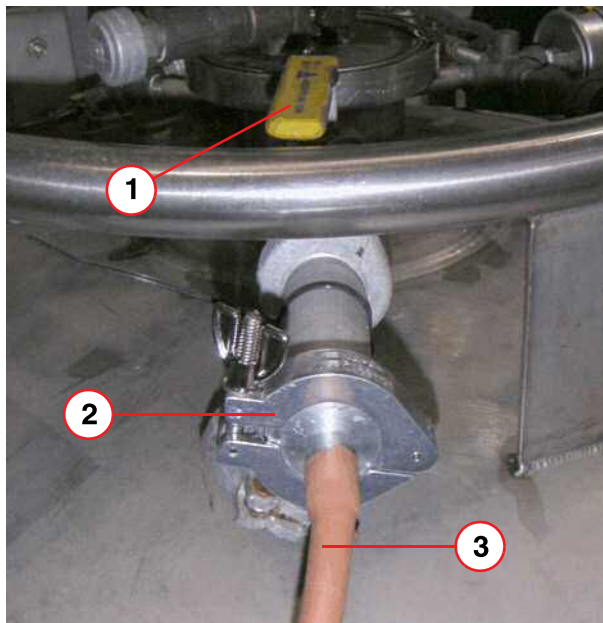


(1) Primary ball valve  
(2) NW25 clamp and seal  
(3) Gas supply line

- Ensure that the **gas supply line** (pos 3) is long enough to allow the dewar to be moved close to the magnet if the magnet is at field.
- **Purge** the helium **gas supply line** of air.
- Connect it to the **primary ball valve** (1) on the dewar, using an **NW25** clamp and seal (pos 2).

- **Open the primary ball valve** (→ 5/ Fig. 10 Page 16) and pressurize the dewar to **0.27bar (4 psi (G))**.

2.2.5 Assemble the service syphon

Fig. 15: Magnet service syphon assembly - top filling

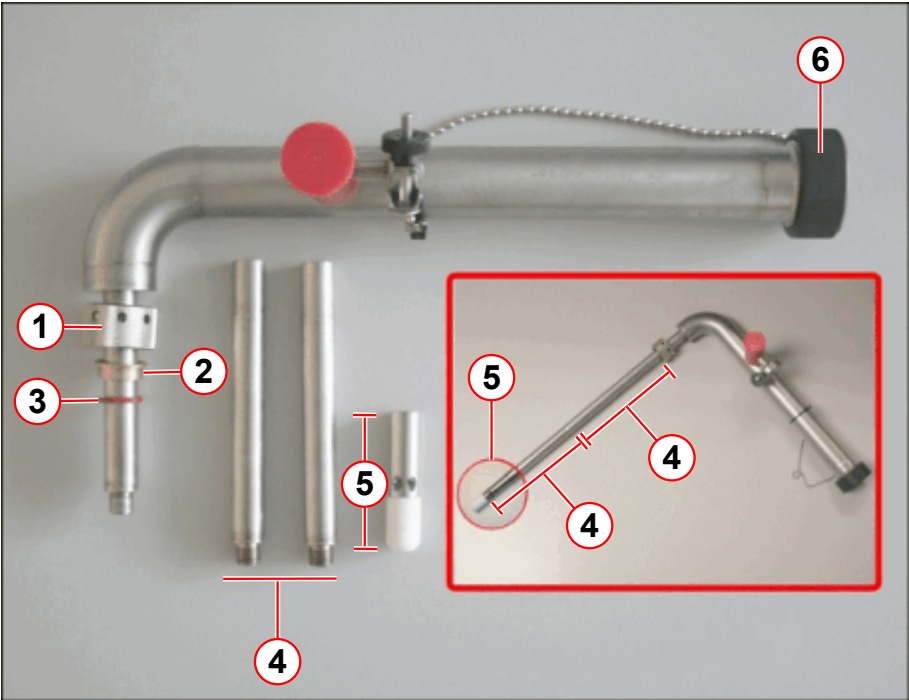


- (1) Ring nut
- (2) Brass flange
- (3) O-ring seal
- (4) 140 mm syphon legs (2 off)
- (5) De-ice and helium fill diffuser and PTFE tip
- (6) Syphon end cap

Tab. 1 OR98 / OR99 / OR103 TOP filling syphon extension legs

| Magnet type                                                                                                                                                          | 140 mm leg | 60 mm De-ice & helium fill diffuser & PTFE tip |
|----------------------------------------------------------------------------------------------------------------------------------------------------------------------|------------|------------------------------------------------|
| Aera (OR98)<br>Skyra (OR99/OR103)                                                                                                                                    | 2          | 1                                              |
| The <b>syphon diffuser</b> should optimally be 50-80 mm above the syphon cone for <b>top filling</b> .<br>If damaged, the <b>PTFE diffuser tip</b> must be replaced. |            |                                                |

- Check the **ceiling height**. (→ Low ceiling height / Page 41)
  - Assemble the **syphon** (only if it can be fitted), or
  - Assemble the **diffuser**, bottom **extension** and support **clamp**.

Fig. 16: Assembling the Ring Nut Seal

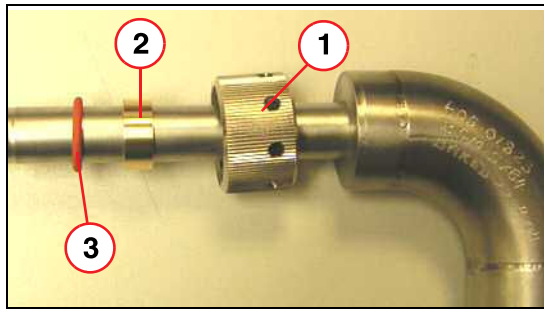

- (1) Ring nut
- (2) Brass spacer
- (3) Silicon O-ring

- Fit the ring nut, brass spacer, and the O-ring to the **service syphon leg**.

Refer to table above (→ Tab. 1 Page 20)

- Connect the **extension legs** and diffuser:
  - Screw the (2 off) **140 mm legs** together.
  - Connect the leg to the **syphon**.
  - Fit the 60 mm **diffuser** and PTFE tip.

## 2.2.6 Fit the magnet service syphon



### CAUTION

Cold head failure.

The cold head seals can fail, and leak to atmosphere.

- » Always ensure the cold head bypass valve is **CLOSED** during helium filling.



### CAUTION

Pressurized helium vessel.


Failure to observe the following safety information may result in physical injury.

- » Do not disconnect any of the fittings without first depressurizing the helium vessel.

### 2.2.6.1 Deactivate the pressure heater

- Access the **Service Software** at the **MRC (MRAWP) Host**.
- Select the menu item **Magnet & Cooling - Magnet Status**.
- At **MSUP Control - Sets the pressure heater mode** select the **OFF** button (→ 1/Fig. 17 Page 22).
- Select **SAVE** (→ 3/Fig. 17 Page 22).

Fig. 17: Service screen Magnet &amp; Cooling - Magnet status - MSUP Control

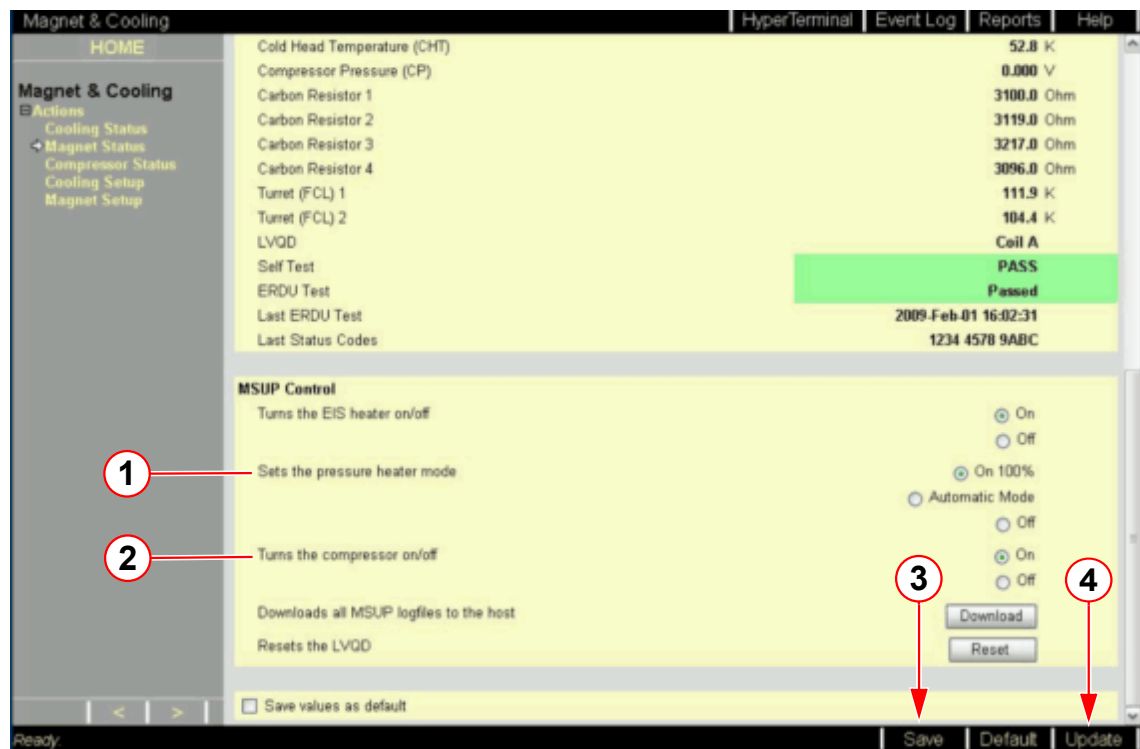

- (1) Sets the pressure heater mode
- (2) Sets the the compressor on/off
- (3) Save button
- (4) Update button



The pressure heater will reactivate within 60 min.

#### 2.2.6.2 Depressurize the helium vessel



#### CAUTION


**Risk of ice formation.**

The magnet bypass valve is only **OPENED** for depressurization of the helium vessel. The bypass valve is **CLOSED** during normal operation.

- » To avoid air ingress, **CLOSE** the magnet bypass valve immediately after the helium vessel is depressurized.




Fig. 18: Bypass valve positions for depressurization

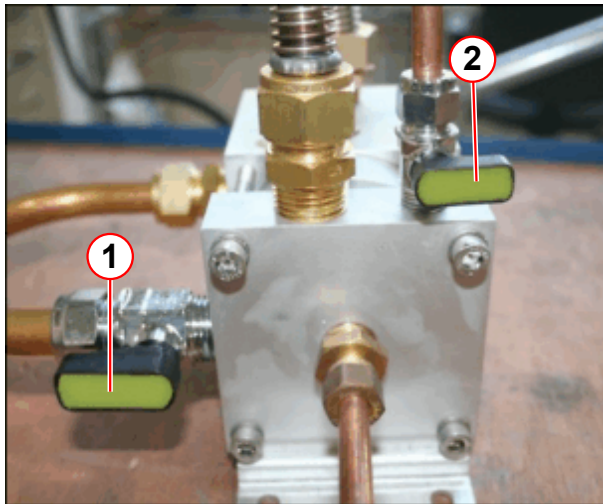

- (1) Magnet bypass valve OPEN  
(2) Cold head bypass valve CLOSED

- **Fully open** the **magnet bypass valve** to **depressurize** the helium vessel.
- Wait until **zero psi** is displayed on the pressure gauge.
- The **pressure switch** will activate.
  - » the **cold head** will turn **off**.
- **Close** the magnet bypass valve to prevent **air ingress**.
- **Warm up** the vent valve assembly.
  - » **Take care** not to damage any components when using a **heat gun**.

Fig. 19: MSUP - Magnet pressure is low

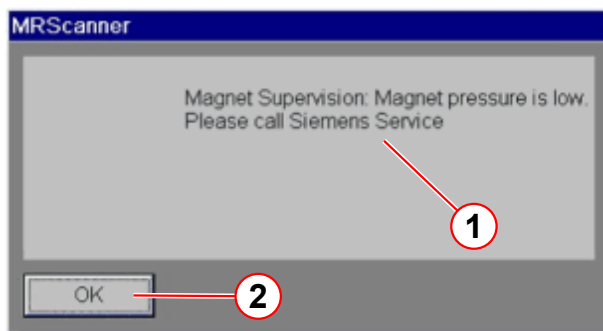


- (1) Pop-up message  
(2) Click OK to clear

On the **MR Scanner**:

- A pop-up **alarm message** is displayed.
    - » **Magnet pressure is low.**
    - » Click **OK** to clear the alarm.
  - The **alarm box** will sound when **pressure falls to alarm limits**.
    - » The yellow **WARNING LED** lights on the alarm box.
    - » The alarm **buzzer** is **activated**
    - » Press **AUDIO ALARM OFF** to silence the alarm.
- Check the **pressure** is positive when the **compressor** switches **OFF**.


Fig. 20: WIKKA Pressure switch - with a small positive pressure

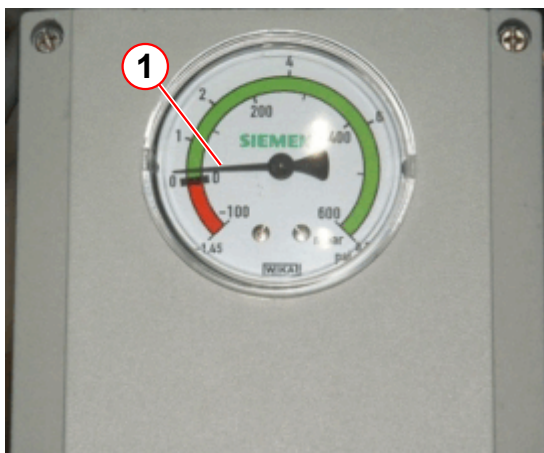

- **Check** the helium vessel has a **small positive pressure**.


- Warm and dry the **syphon blanking plug** using a heat gun to ease removal.  
(→ 1/Fig. 22 Page 24)

Fig. 21: Syphon end cap

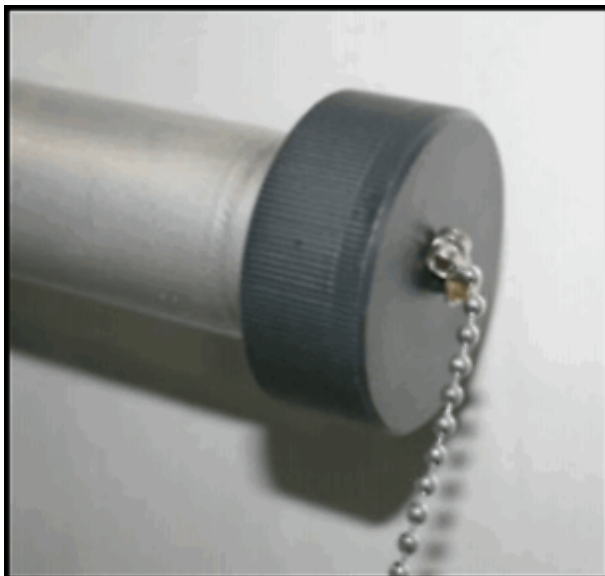

- Remove the syphon **end cap**.

Fig. 22: Fit/remove the syphon blanking plug

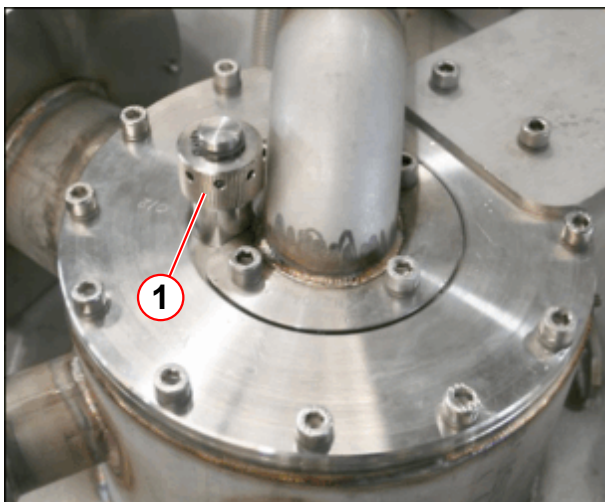

- Remove the **syphon blanking plug** using the syphon 'C' spanner.



### CAUTION


Risk of cold burns.

Cold gas will exhaust from the syphon port on opening.

- » Stand clear when opening.
- » Wear protective gloves and safety glasses.



Fig. 23: Helium syphon

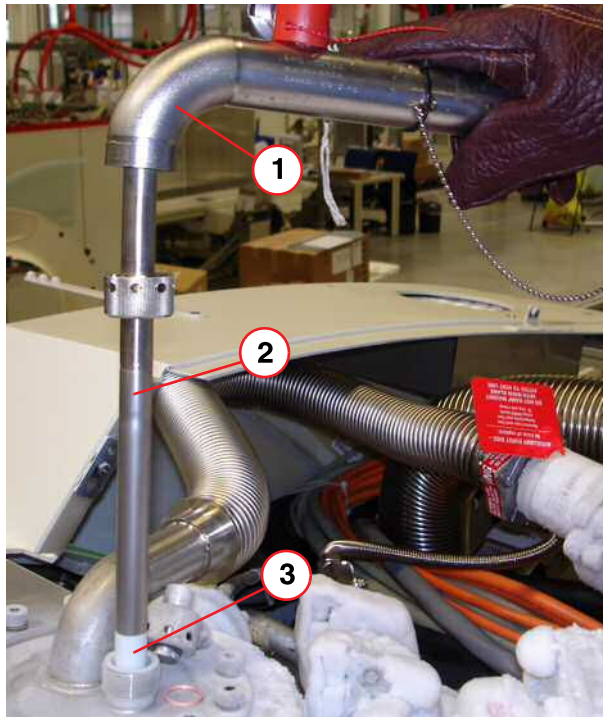

- (1) Service Syphon  
 (2) Two Extension Legs  
 (3) PTFE Syphon Diffuser

- **Insert** the service syphon into the port without delay.
- Slowly lower the service syphon fully into the magnet.
- Tighten the **syphon port ring nut** by hand, then apply a half turn using the syphon C-spanner to form a **gas tight seal**.



There is a risk of quenching the magnet if the syphon is inserted too quickly.

- **Re-fit** the syphon **end cap**. (→ Fig. 21 Page 24)



If a helium fill is showing poor fill efficiency, or no increase on the helium level meter, the transfer must be stopped and the syphon port inspected for ICE.

### 2.2.7 Precool the service syphon



The helium dewar must be positioned close to the magnet with the dewar syphon fitted, ready for use, (Refer to (→ Dewar preparation / Page 15)).



If the temperature of the magnet is above 4.2K, a dense cascading plume will not form on cooling the service syphon.

**CAUTION****Risk of quench**

Warm helium gas can cause a quench.

- » If the magnet is at field, it will quench if both syphons are not pre-cooled correctly.
- » Ensure syphons are correctly pre-cooled before joining together.
- » Frost or condensation on the body of either syphon is an indication of poor vacuum. Stop the transfer immediately - the syphon may need replacing.
- » Frost on the syphon coupling is acceptable.

Cool down the **service syphon** as follows:-


Fig. 24: Syphon end cap



- Ensure the **syphon end cap** is fitted and the **bypass valve** is closed.

- **At the MRC, reactivate** the **pressure heater** to achieve a positive pressure of about **0.3 psi (G)** in the magnet. (→ Reactivate the pressure heater / Page 34)
- **Deactivate** the pressure heater (→ Deactivate the pressure heater / Page 21).

Fig. 25: Pre-cooling the syphons

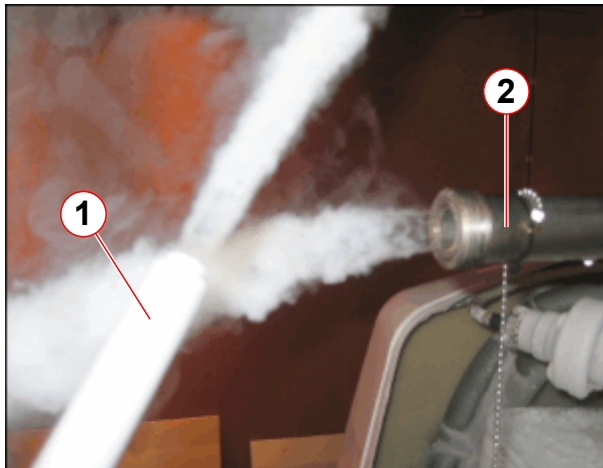

- (1) Dewar syphon  
(2) Service syphon

- Remove the **syphon end cap**.
- Allow the syphon to vent until **fully cold**.
  - a “cascading plume” of helium gas appears.

**CAUTION**

Risk of cold burns.

Cold gas will exhaust from the syphon port on opening.

- » Stand clear when opening.
- » Wear protective gloves and safety glasses.


### 2.2.8 Liquid helium transfer - top fill



Before helium transfer, refer to (→ Helium pressures at altitude / Page 42)

- Check the **dewar** gas supply pressure is set at **0.27bar (4psi (G))**.
- As the **service syphon** nears the point of being fully cold:
  - open the **dewar syphon flow valve** half way (2 turns). (→ Fig. 27 Page 29)
  - **precool the syphon** again until a dense white “cascading plume” appears.

Fig. 26: Dewar and service syphon coupling

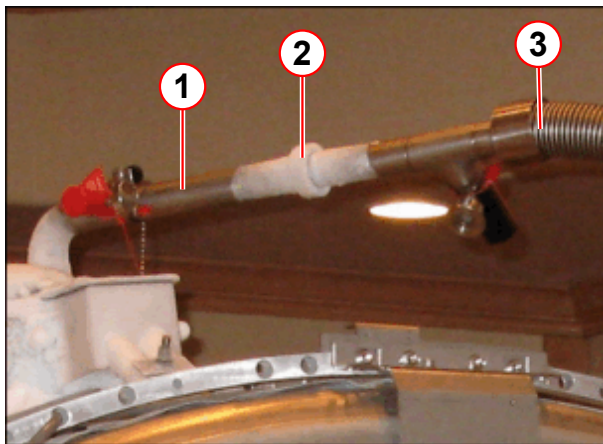

- (1) Service syphon
- (2) Lock nut
- (3) Dewar syphon

- When **both** syphons are giving a dense white "cascading plume" of helium gas, connect them together and secure the lock nut.
- OPEN the magnet bypass valve. This will:
  - » Improve efficiency.
  - » Improve the filling speed.




## CAUTION

### Risk of quench

Warm helium gas can cause a quench.

- » If the magnet is at field, it will quench if both syphons are not pre-cooled correctly.
- » Ensure syphons are correctly pre-cooled before joining together.
- » Frost or condensation on the body of either syphon is an indication of poor vacuum. Stop the transfer immediately - the syphon may need replacing.
- » Frost on the syphon coupling is acceptable.

Fig. 27: Dewar syphon flow valves - open and closed positions

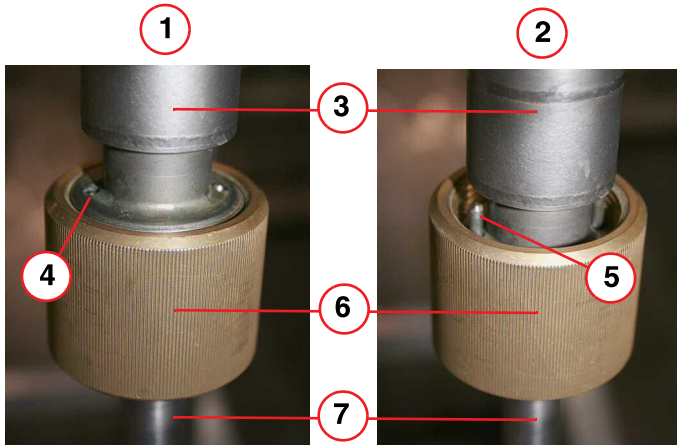

- (1) FULLY OPEN position
- (2) FULLY CLOSED position
- (3) Stainless steel body
- (4) Guide pin flush with body
- (5) Guide pin
- (6) Brass nut
- (7) Syphon leg

Fig. 28: Helium filling



- Fully opening the **dewar syphon flow valve**.
  - » The flow valve is **fully open** when the two **guide pins** at the top of the brass nut are flush with the top of the stainless steel body (→ 4/ Fig. 27 Page 29).
  - » From fully closed to fully open, is **4 turns**.
- Ensure the transfer dewar pressure is **4 psi (A)** at normal altitude. Refer to (→ Helium pressures at altitude / Page 42).

### 2.2.8.1 Perform a helium level check



During installation, if the MRC is not available:


1. Take a helium sample using the MSUP. (→ Helium level check at the MSUP / Page 31)

### Helium level check at the MRC

- Check the liquid helium level by using Service Software on MRC (MRAWP) Host.
  - Check the liquid helium level using the MRC.
  - Access the Service Software at the MRC.
  - Select the menu item Magnet & Cooling - Magnet Status.

- Press the Start button (→ 3/ Fig. 29 Page 30) to initiate the helium measurement.  
(→ 1/ Fig. 29 Page 30)

Fig. 29: Checking liquid helium and magnet pressure levels

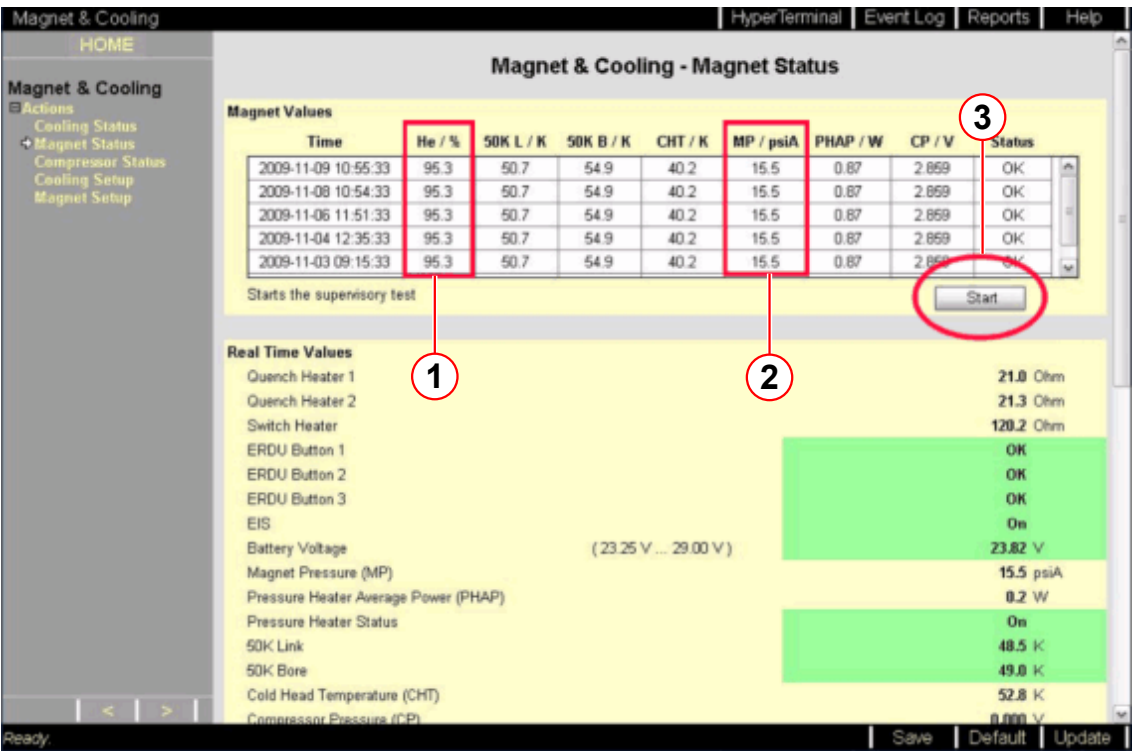


- (1) Liquid helium level
- (2) Magnet pressure
- (3) Start button

**i** Data is only accurate from the last refresh of the page, or from when the page was opened last.


## Helium level check at the MSUP

Fig. 30: MSUP 2306 - depress SAMPLE to activate

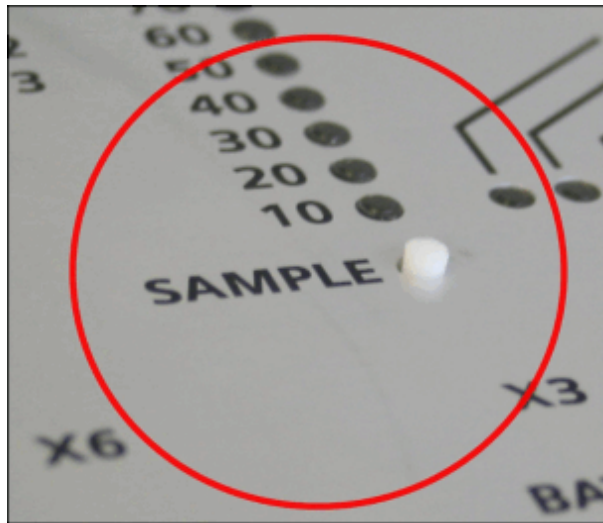

- Locate the SAMPLE switch on the MSUP lid.
- Depress the sample switch for 2 s, and release.

Fig. 31: MSUP - Helium level sample

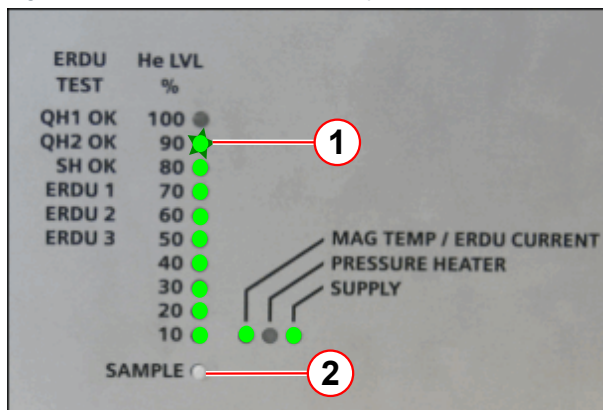

- (1) Sample switch  
(2) Helium level

- The green LEDs illuminate and display the actual liquid helium level percentage.
- If the 90% LED is flashing, STOP the transfer. Refer to (→ Tab. 2 Page 31).

Tab. 2 MSUP Helium display

| MSUP LEDs                      | Helium level | Comments or actions            |
|--------------------------------|--------------|--------------------------------|
| 0 - 80% ON                     | As displayed | Continue to fill as required.  |
| 90% ON                         | As displayed | Stop transfer.                 |
| *90% FLASHING<br>10 - 80% ON   | 95% - 97%    | STOP the transfer IMMEDIATELY! |
| *100% FLASHING<br>10% - 90% ON | 97% and over | Magnet is OVERFILLED.          |





\* The LEDs return to ON after approximately 30 seconds.



To avoid overfilling, the helium transfer must be stopped as soon as the 90% LED illuminates.

- Repeat the measurement every 5 minutes. An increase in level should be seen within 20 minutes.
  - » If not, stop the transfer and investigate the reason for poor transfer efficiency.



Opening the bypass valve will help to reduce the back pressure, enabling more efficient helium transfer.



Maintain a pressure of about 0.27 bar (4 psi (G)) at the helium dewar. Too much or too little pressure will result in inefficiency.

At 4 psi (G), the transfer rate is approximately 10 liters per minute. **Do not leave the transfer unattended.** A 250 liter dewar will empty in 25 - 30 minutes.

If frost or condensation forms on the body of either syphon, it is an indication of poor vacuum. **Stop the transfer immediately - the syphon may need to be replaced.**

## NOTICE

**Risk of quench.**

**Warm helium gas will cause a quench.**

- » **Stop the transfer immediately in case you hear whistling (a high pitched sound) from the syphon, or the dewar pressure falls.**

### 2.2.9 Removing the service syphon

- **Depressurize** the helium vessel. (→ Depressurize the helium vessel / Page 22)


Fig. 32: Warm up the syphon joint with heat gun.



- **Warm** the connecting **syphon joint** to room temperature (using a **heat gun**) to remove ice.

- **Warm** the syphon **entry valve** to room temperature and dry.

Fig. 33: Syphon end cap

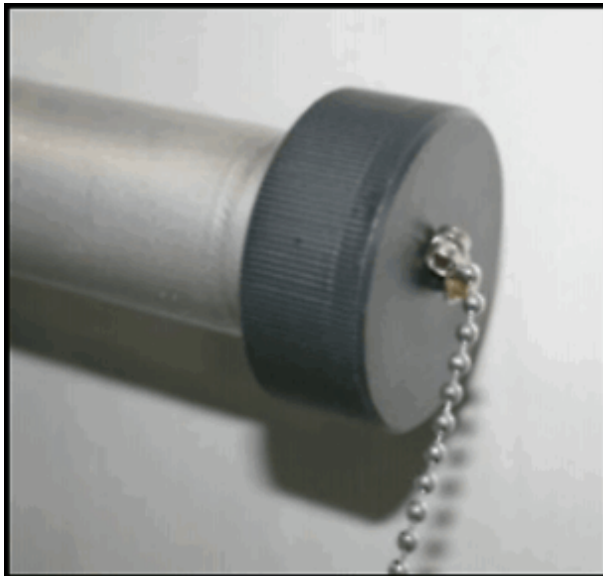

- **Warm** the syphon **end cap** ready for fitment onto the syphon.

- **Disconnect** the syphons.



#### CAUTION

Risk of cold burns.


Cold gas will exhaust from the syphon port on opening.

- » Stand clear when opening.
- » Wear protective gloves and safety glasses.

- **Close off** the end of the service syphon with the **syphon end cap**. (→ Fig. 33 Page 33)


- **Immediately remove** the service syphon from the magnet.

Fig. 34: Syphon blanking plug

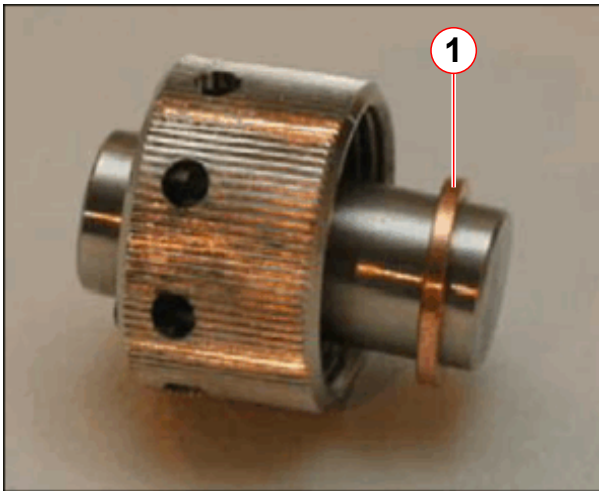

(1) Copper seal

Fig. 35: Fit/remove the syphon blanking plug

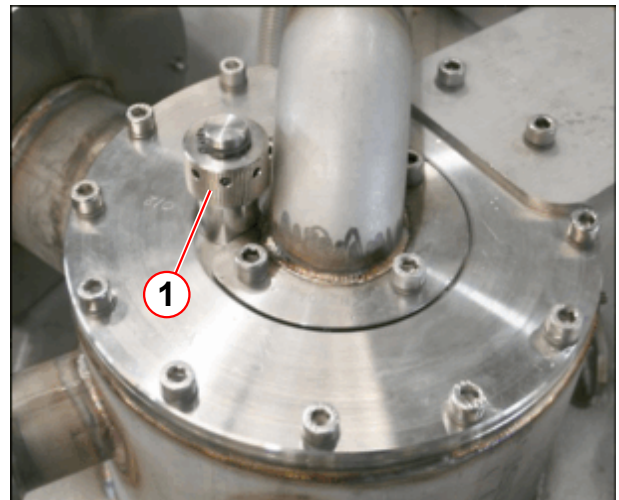

- **Immediately fit the syphon blanking plug** using a new copper seal.  
(→ Fig. 34 Page 34)
- **Warm up the threads** on the plug and **carefully** tighten with the syphon 'C' spanner.
  - There might be ice on the threads - **do not use undue force**. Damaged threads will result in a leak.
- Warm the **top plate** and **syphon plug** using a heat gun. **Re-tighten** the plug.



### CAUTION

**Danger of air ingress.**

**Internal ice blockage could cause excessive pressure.**

- » Warm up syphon port plug sufficiently to ensure there is no ice on the threads when tightening.

- **Pressurize** the helium vessel to **15.5 psi (A)**.

#### 2.2.9.1 Reactivate the pressure heater

- Access the **Service Software** at the **MRC (MRAWP) Host**.
- Select the menu item **Magnet & Cooling - Magnet Status**.
- At **MSUP Control - Sets the pressure heater mode** select the **On 100%** button.  
(→ 1/Fig. 17 Page 22)
- Select **SAVE**. (→ 3/Fig. 17 Page 22)

- **At Real-Time Values** - Monitor the **magnet pressure**.
  - Select **Update** to refresh data until pressure reaches  $\geq 15.5$  psi (A).




Data is only accurate from the last refresh of the page, or from when the page was opened last.



During pressurization, the cold head will switch on automatically.  
If the cold head was turned off at MRC, select (→ 2/Fig. 17 Page 22)

Fig. 36: Check magnet pressure level

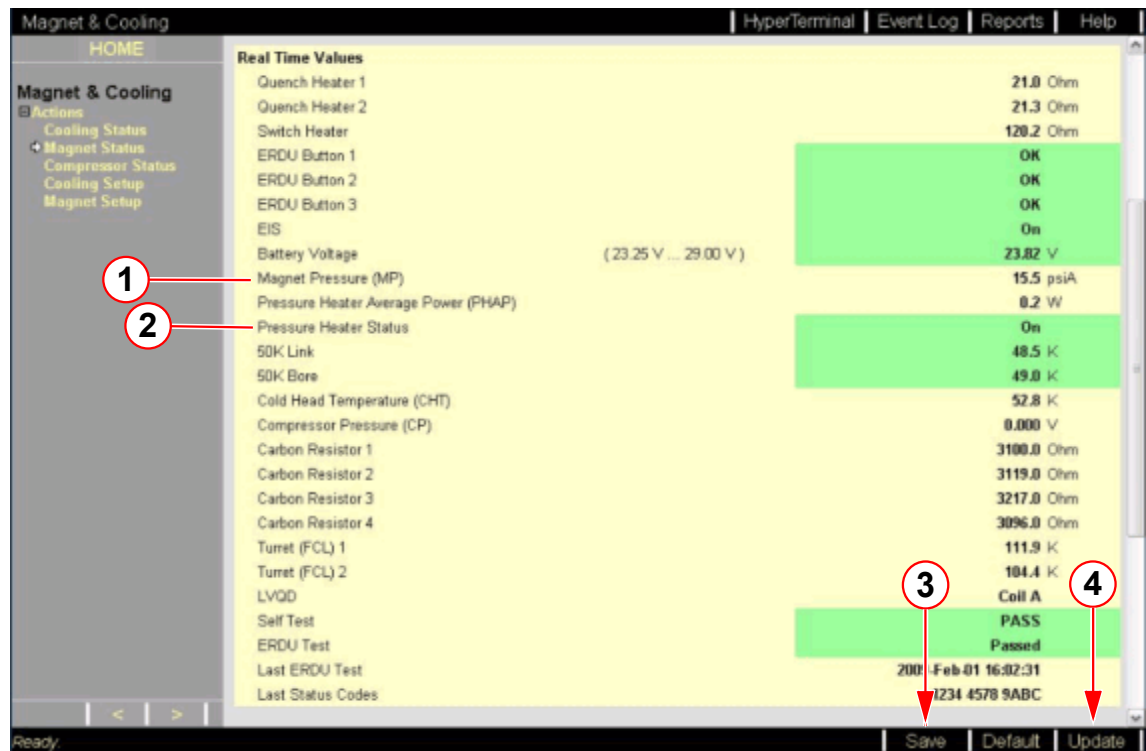

- (1) Magnet pressure (MP)
- (2) Pressure heater status
- (3) Save button
- (4) Update button

- **At MSUP Control** - Sets the pressure heater mode select the **Automatic** mode button. (→ 1/Fig. 17 Page 22)
- Select **SAVE**. (→ 3/Fig. 36 Page 35)

### 2.2.9.2 Leak check

- After helium fill has been completed, **leak check** the **turret**.



Wait at least 2 hours after completion of a helium fill before attempting to ramp the magnet.

- **Move the dewar** out of the imaging suite into a **well ventilated area**.

- Complete **depressurizing** the dewar by opening the **primary ball valve**. (→ 5/Fig. 10 Page 16)
- When at the gas pressure is at **zero**, **close** the **primary ball valve** on the dewar.
- **Remove** the dewar syphon.
- **Immediately** close the **syphon entry valve** on top of the dewar. (→ 4/Fig. 10 Page 16)
- **Open** the **low-pressure dewar relief valve**. (→ 1/Fig. 10 Page 16)
- **Warm** both **syphons** with the heat gun **until dry**, and store safely.

### 2.2.10 Changing supply dewars during transfer




Whistling of the syphon (a high pitched sound), or a noticeable drop in dewar pressure are indications of an empty dewar.

If the supply dewar empties before the magnet is filled, it will be necessary to replace the supply dewar with a full one.

Ensure that a full dewar is available and depressurized, ready to receive the syphon.

Fig. 37: Dewar syphon fitment

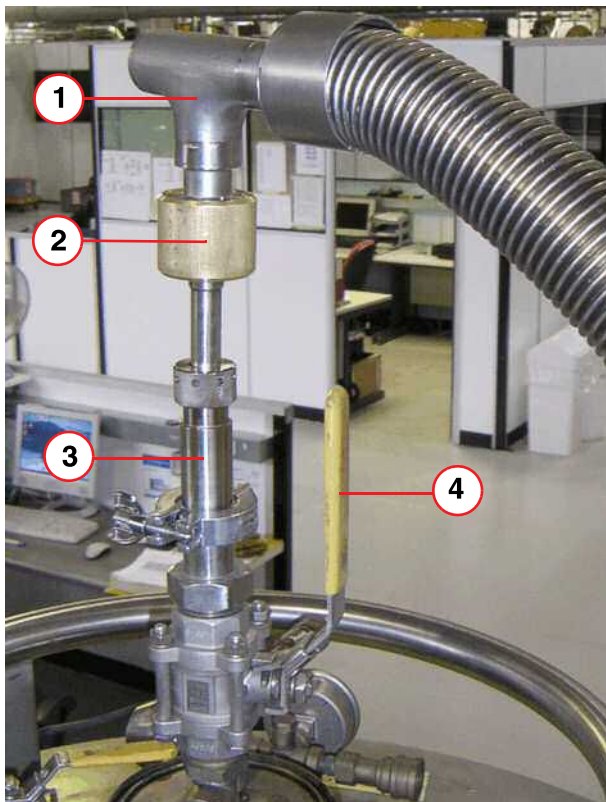

- When the dewar empties, stop the transfer **immediately** by closing the **dewar syphon flow valve** (pos 2).

- (1) Syphon
- (2) Syphon flow valve
- (3) Dewar syphon adaptor
- (4) Syphon entry valve

**CAUTION**


**Risk of quench.**

Failure to close the valve immediately, when the magnet is at field, will result in a quench.

» Close the valve immediately.

- Turn **off** the dewar gas supply:
  - remove the **gas connection**.
  - vent the dewar down to atmospheric pressure.
  - **Ensure the area is well ventilated.**
- **Close the primary dewar ball valve.** (→ 5/ Fig. 38 Page 37)

Fig. 38: Dewar Ports (typical)

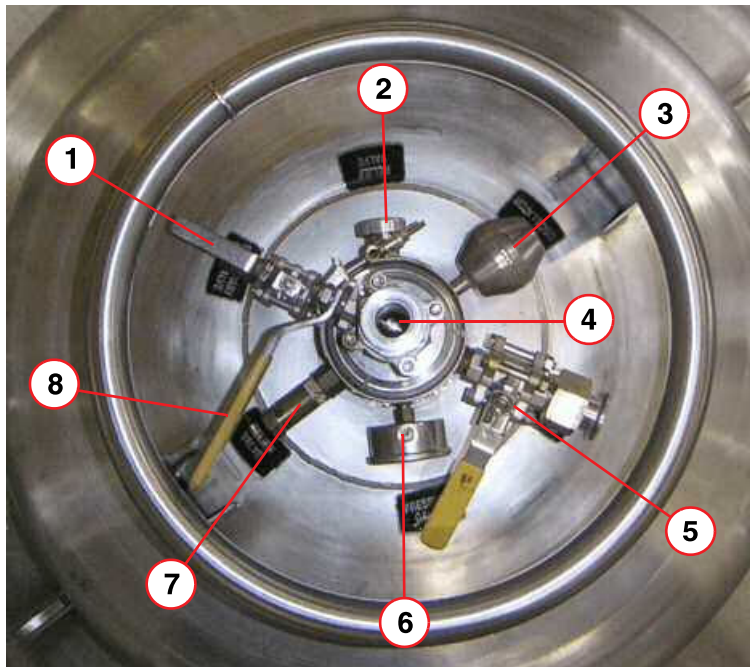

- (1) Low pressure relief valve (shown open)
- (2) Relief valve
- (3) Anti-oscillator
- (4) Syphon entry port (shown closed)
- (5) Primary ball valve (shown closed)
- (6) Pressure gauge
- (7) Relief valve
- (8) Syphon entry valve

- Using a heat gun, **warm the syphon connecting joint** to room temperature.



- Disconnect the syphons.

**CAUTION**

Risk of cold burns.

A "cascading plume" of cold gas will exhaust from the syphons when they are disconnected.

- » Stand clear of the syphon joint when disconnecting.
- » Wear protective gloves and safety glasses.

- Remove any moisture from the end of the **service syphon** with paper towelling.
- Fit the **service syphon end cap**. (→ Fig. 33 Page 33)

**CAUTION**

Danger of air ingress.

Internal ice blockage could cause excessive pressure.


- » Ensure the syphons are capped and plugged immediately and securely.

- Move the empty dewar close to the full dewar.
- Vent both storage dewars down to **atmospheric pressure**.
  - » This must be done in a well ventilated area - NOT in the imaging suite.
- Open the **low pressure relief valve** on the **empty** dewar. (→ 1/ Fig. 38 Page 37)
- Remove the dewar **syphon** and **adaptor** from the empty dewar and **immediately** insert the syphon leg into the **full dewar**.
- Close the **syphon entry valve** (→ 4/ Fig. 38 Page 37), of the **empty** dewar. and **open the low pressure relief valve**. (→ 1/ Fig. 38 Page 37)
- Close the **low pressure relief valve** of the **full** dewar.
- Purge the gas supply line.
- Connect to the full dewar and **open** the dewar **primary ball valve**.
- Set the **gas pressure** to 0.27bar (4psi (G)).
- Precool the **dewar syphon**. (→ Precool the dewar syphon / Page 18)
- Move the dewar next to the magnet system.



- **Precool the service syphon.**

Fig. 39: Dewar and service syphon coupling

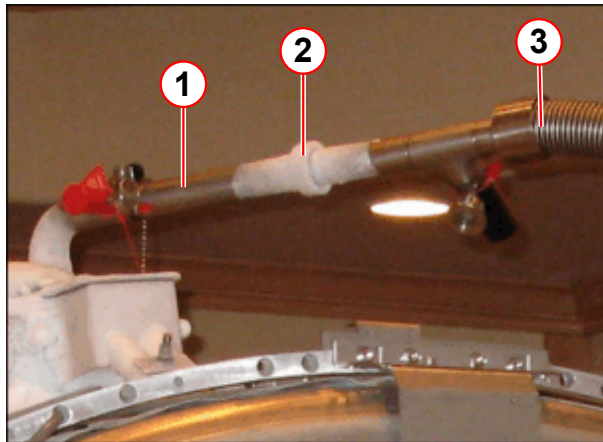

- (1) Service syphon
- (2) Lock nut
- (3) Dewar syphon

- **Continue with the transfer.**

- When **both** syphons are giving a dense white "cascading plume" of helium gas, **connect** them together and **secure** the lock nut.

### 2.2.11 Stopping the helium transfer

Fig. 40: Dewar manometer



- As soon as the helium level reads the required level or at a maximum **95%**, or the dewar is empty, **stop the transfer** by closing the flow valve at the dewar syphon. **Do not fill above 95%.** (→ Helium level check at the MRC / Page 29)
  - Whistling (high pitched sound) of the syphon or noticeable drop in dewar pressure are indications of an empty dewar.
- The **dewar syphon flow valve** must be closed **immediately** the dewar empties - **if left open there is a risk of quenching.** (→ 6/ Fig. 27 Page 29)
- Set the dewar gas pressure regulator to **zero**.
- Close the dewar **primary ball valve.** (→ 1/ Fig. 14 Page 19)
- **Remove the gas line from the dewar.**
- Refer to (→ Changing supply dewars during transfer / Page 36) if the magnet is not full.

## NOTICE

**Risk of quench.**

Warm helium gas will cause a quench.


- » If the magnet is still at field, stop the transfer immediately in case you hear whistling (a high pitched sound) from the syphon, or the dewar pressure drops.

### 2.2.12 Final work steps

After transfer, check:-


- Press the **LVQD reset button** to ensure the LVQD shows **waiting trigger**.

Fig. 41: LVQD reset button

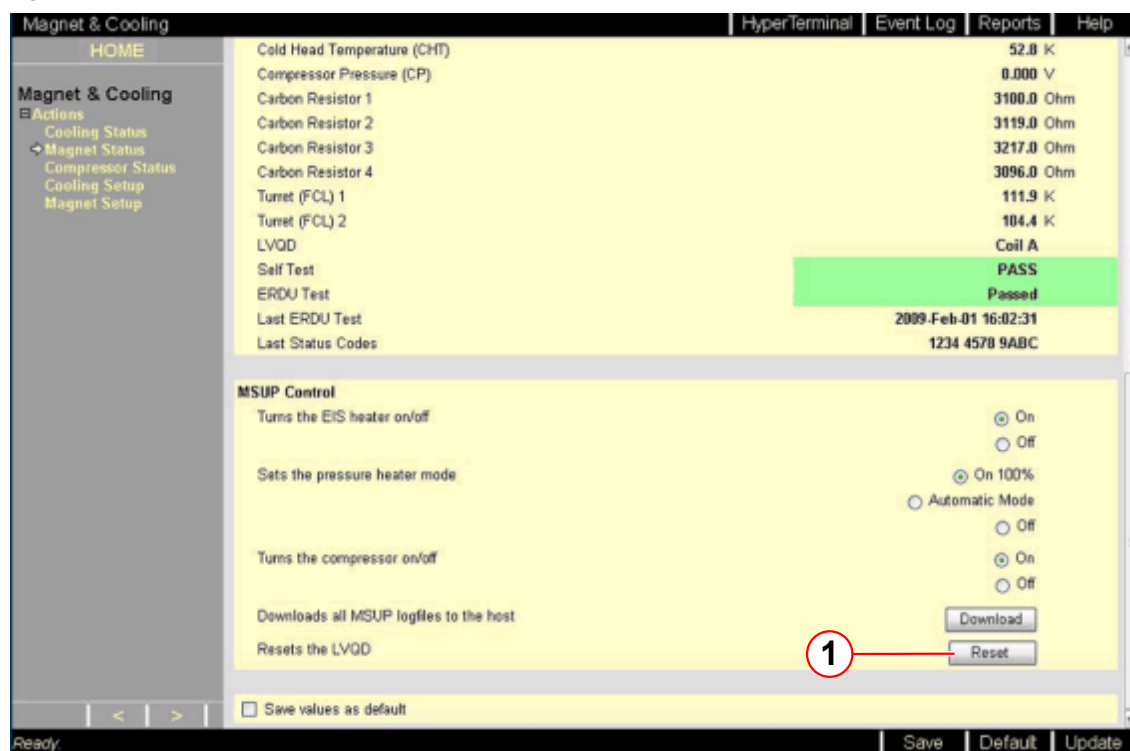

- The **magnet bypass valve** is closed.
- The **pressure heater** is operating correctly.

Fig. 42: MSUP - Pressure heater LED

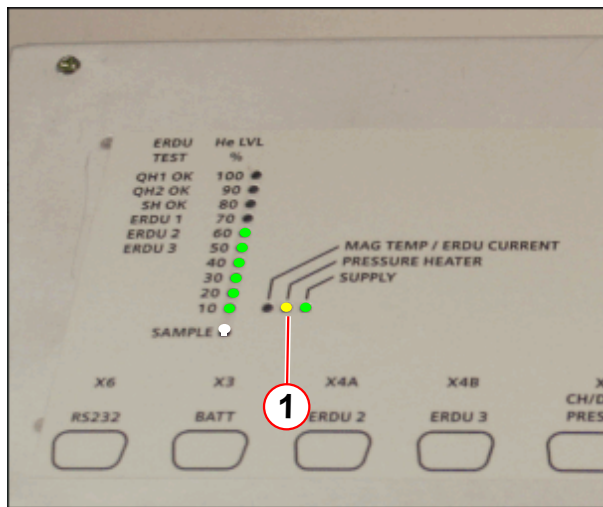

(1) Pressure heater LED - yellow

- The pressure heater **LED** on the **MSUP** is illuminated, or switching on and off.


- **No alarm** signals are evident at the **MSUP** or the **alarm box**.
- The **compressor** is switched **on** and **running**.
- The **magnet pressure** is stable at **15.5psi (A)**.
- The **turret** and **venting** system are **warm** and **dry**.
- The filling **port** and **venting** system have **no leaks**. (Refer to **Venting Leak Check procedure**).
- The **4K cold head flange** has **no leaks**.
- The used **dewars** are correctly **closed**, with their low pressure relief valves **open**.

### 2.2.13 Low ceiling height

Where magnets have been installed into rooms with a **low ceiling height**, it can be difficult to fit the assembled service syphon. To avoid dropping the extension leg(s) down into the syphon port, the **syphon leg support clamp** must be used (→ 3/ Fig. 43 Page 42).

- Check the **ceiling height**:
  - Assemble the syphon (only if it can be fitted), or
  - Assemble the diffuser, bottom extension and support clamp

Fig. 43: Syphon configuration for low ceiling height

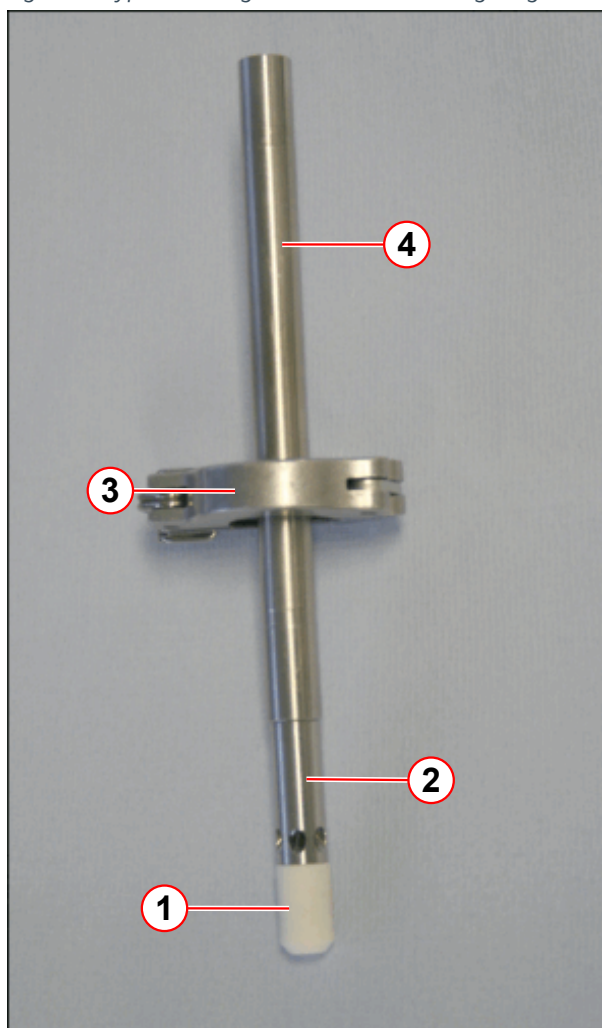

- (1) PTFE extension tip (top fill)
- (2) Helium fill diffuser tip (top fill)
- (3) Syphon leg support clamp
- (4) Extension leg

#### Refer to Magnet Service syphon kit

- Remove the **syphon end cap**.
- Ensure the **ring nut, brass spacer and O-ring** are fitted correctly.
- Fit the **diffuser** to the extension leg.
- Fit **support clamp** to the extension leg.
- **Warm** the syphon **blanking plug** to ease removal.
- **Remove** the syphon **blanking plug** using the C-spanner.
- **Lower** the assembly into the **neck**.
- **Screw** on the second **extension leg**.
- **Move** the syphon **clamp** up the extension leg and lower extensions into the entry port, to allow room for the syphon body.
- **Screw** the **diffuser leg** to the syphon **body**.
- **Remove** the syphon **clamp** and slowly **lower the syphon fully** into the syphon entry.
- **Screw** the **syphon port ring** nut into place and **tighten** to make a gas-tight **seal**.
- **Re-fit** the syphon **end cap**.



Use the support clamp for removing the syphon.

## 2.2.14 Helium pressures at altitude

### 2.2.14.1 Helium transfer with magnet bypass valve closed

As a function of altitude, the graph(s) show the **gauge pressure** in the magnet.


The optimum pressure differential between magnet and dewar for a helium transfer is 4 psi. To calculate this:

- Determine the **magnet's altitude**

- Determine the **magnet gauge pressure** from the graphs below.

#### 2.2.14.2 OR98 magnet - 16 psi (A) valve

Fig. 44: Effect of a change in altitude on gauge pressure (16 psi (A))

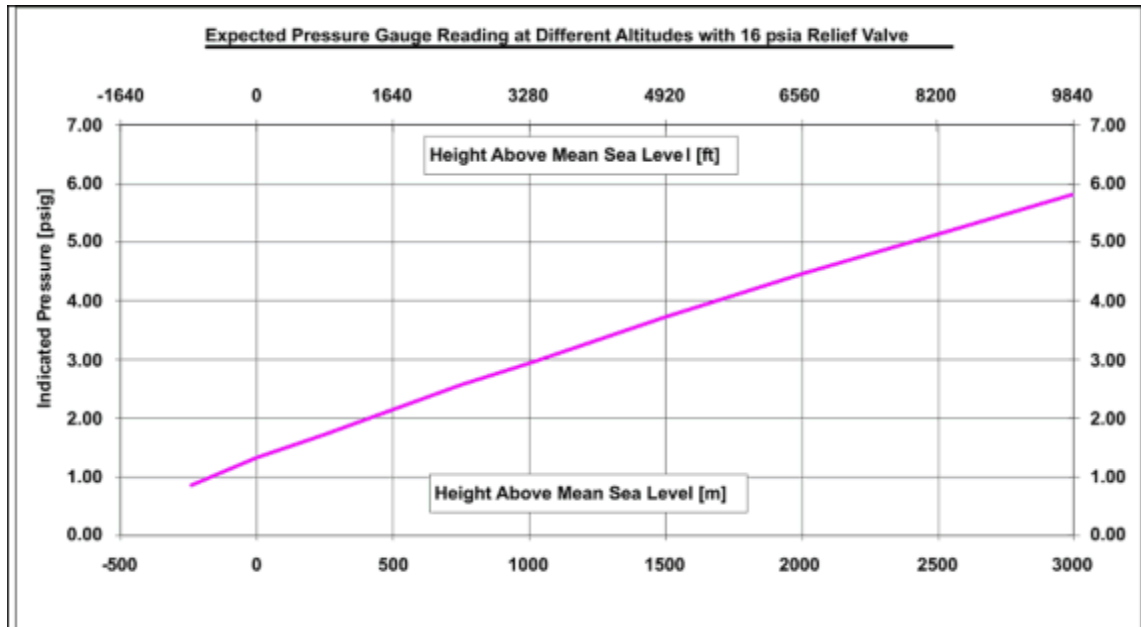

Tab. 3 Dewar transfer pressures (typical)

|                        |   |      |      |      |
|------------------------|---|------|------|------|
| Altitude (M)           | 0 | 1000 | 2000 | 3000 |
| Dewar pressure psi (G) | 5 | 7.0  | 8.5  | 9.5  |

- Add 4 psi to this number to calculate the **required transfer pressure**

» **Example of 16 psi (A) vent valve assembly**


Altitude (1000m) = 3.0 psi (G) + 4 psi = 7.0 psi (G)



Maximum dewar pressure is 10.5 psi (G).

## 2.2.14.3 OR99/OR103 magnet 17 psi (A) valve

Fig. 45: Effect of a change in altitude on gauge pressure (17psi (A))

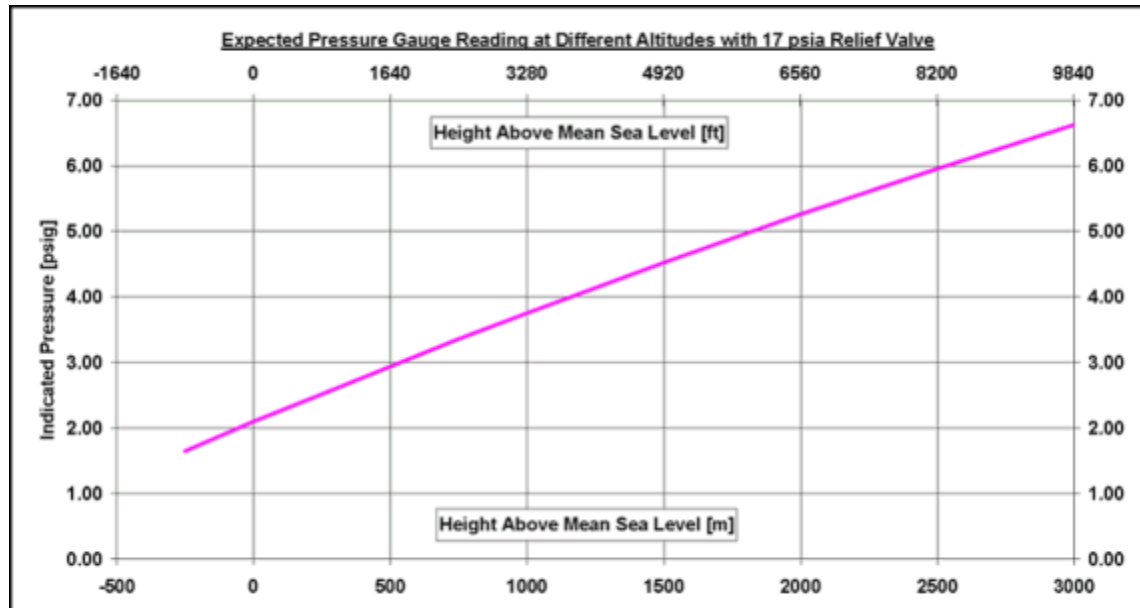

Tab. 4 Dewar transfer pressures (typical)

|                        |   |      |      |      |
|------------------------|---|------|------|------|
| Altitude (M)           | 0 | 1000 | 2000 | 3000 |
| Dewar pressure psi (G) | 6 | 7.75 | 9.25 | 10.5 |

- Add 4 psi to this number to calculate the **required transfer pressure**
  - » Example of 17 psi (A) vent valve assembly  
Altitude (1000m) = 3.75 psi (G) + 4 psi = 7.75 psi (G)



Maximum dewar pressure is 10.5 psi (G).

## 2.2.15 Verifying the electronic low pressure switch

This procedure may be carried out in conjunction with any other process involving depressurization of the helium vessel, i.e. fitting the service syphon, cold head replacement. Pressure control of the helium vessel is managed by the **MSUP**, via the **Pressure heater control unit**. The low pressure switch must be checked **whenever the magnet is depressurized**.

When pressure in the helium vessel drops to **within 0.2 psi (G)** above atmospheric pressure, the **electrical pressure switch** will open and turn **off** the compressor.



Atmospheric pressure may differ, depending on the altitude and local weather conditions.

As the helium vessel is **pressurized**, the switch closes and the **compressor** will **re-start**.



The pressure switch is factory calibrated to operate at 0.2 psi (G) above atmospheric pressure.




Perform this check after helium filling when the magnet syphon is \*removed.

The filling port is closed and the refrigerator switched on.

Note: \*Mobile systems have a fixed syphon.

- Access the **Service Software** at the **MRC (MRAWP) Host**.
- Select the menu item **Magnet & Cooling - Magnet Status**.
- Check **magnet pressure** is normal (**15.5 psi (A)**).

Fig. 46: Check magnet pressure level

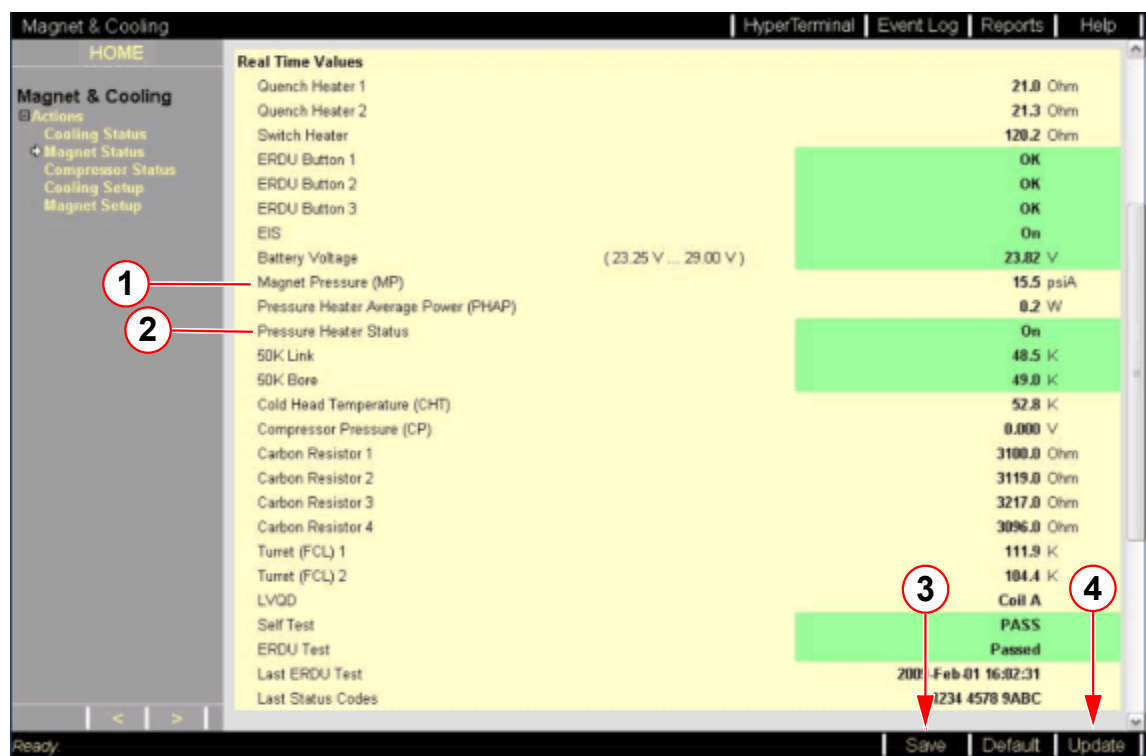

- (1) Magnet pressure (MP)
- (2) Pressure heater status
- (3) Save button
- (4) Update button


- Check the **compressor** is switched **on** and **running**.



### Deactivate the pressure heater

- At **MSUP Control** - Sets the pressure heater mode select the **OFF** button.  
(→ 1/ Fig. 47 Page 46)

Fig. 47: Service screen Magnet & Cooling - Magnet status - MSUP Control

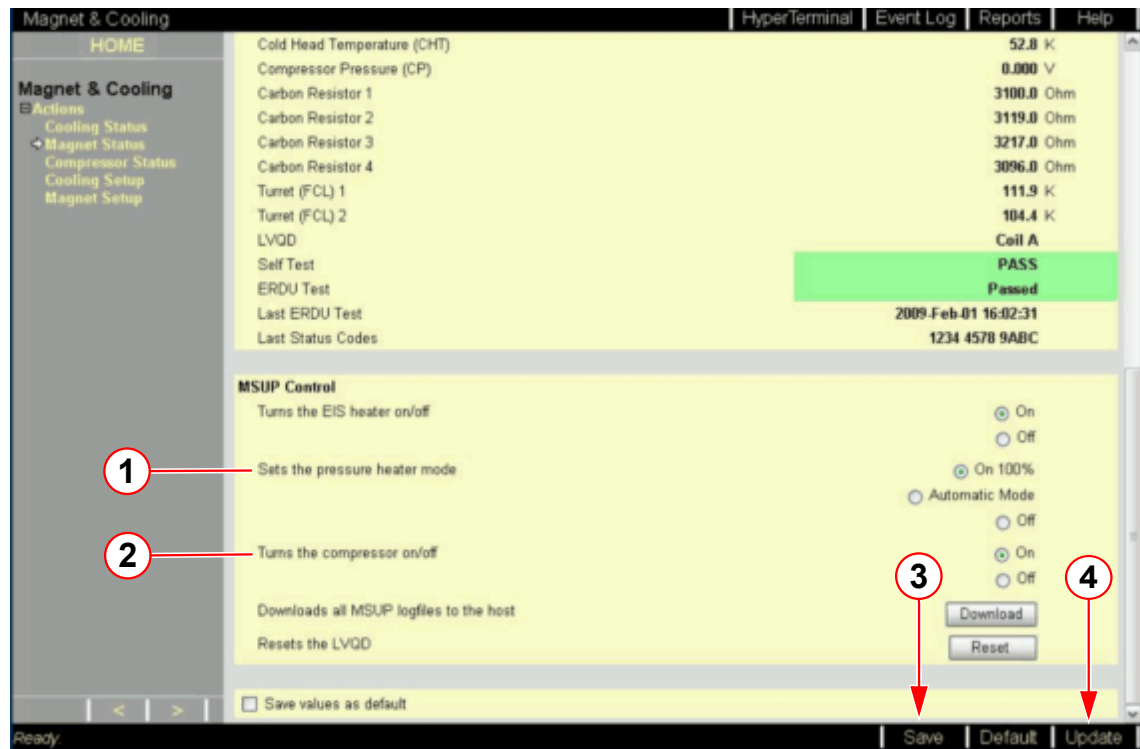


- (1) Sets the pressure heater mode
- (2) Sets the the compressor on/off
- (3) Save button
- (4) Update button

#### 2.2.15.1 Depressurize the helium vessel

- Depressurize** to within **0.2 psi (A)** of atmospheric pressure.


Fig. 48: Bypass valve positions for depressurization

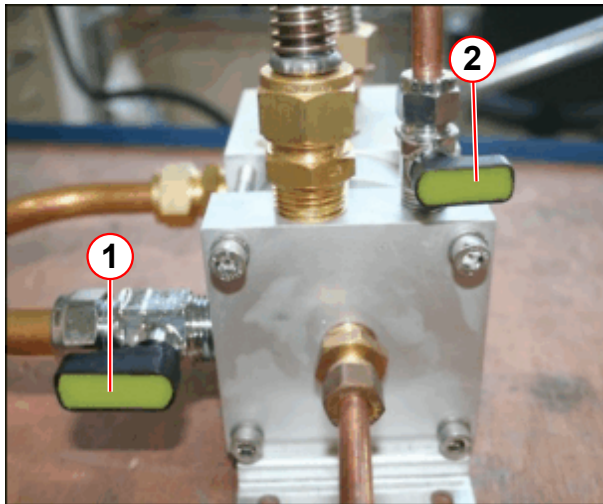

- (1) Magnet bypass valve OPEN  
(2) Cold head bypass valve CLOSED

- **Fully open the magnet bypass valve to depressurize the helium vessel.**
- The **pressure switch** will activate when system pressure is within **0.2 psi (A)** of atmospheric pressure.
  - » this turns the **cold head OFF**.
- **Close the magnet bypass valve to prevent air ingress.**
- **Warm up the vent valve assembly.**
  - » **Take care** not to damage any components when using a **heat gun**.

Fig. 49: MSUP - Magnet pressure is low

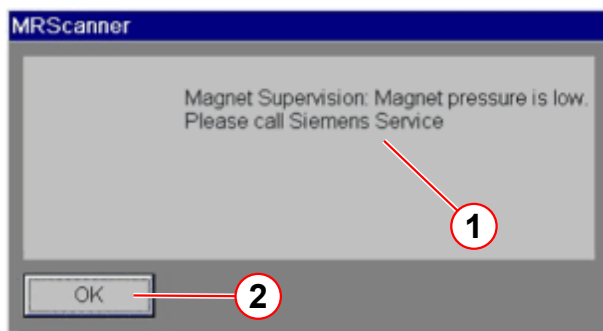

- (1) Pop-up message  
(2) Click OK to clear

On the **MR Scanner**:

- A pop-up **alarm message** is displayed.
    - » **Magnet pressure is low.**
    - » Click **OK** to clear the alarm.
  - The **alarm box** will sound when **pressure falls to alarm limits**.
    - » The yellow **WARNING LED** lights on the **alarm box**.
    - » The alarm **buzzer is activated**
    - » Press **AUDIO ALARM OFF** to silence the alarm.
- Check the **pressure is positive** when the **compressor** switches **OFF**.
  - **Repressurize** the helium vessel.<sup>1</sup>

#### 2.2.15.2 Reactivate the pressure heater

- At **MSUP Control - Sets the pressure heater mode** select the **On 100%** button.  
(→ 1/Fig. 47 Page 46)
- Click **UPDATE**. (→ 4/Fig. 47 Page 46)
- Check that the **compressor** activates **automatically** as the pressure rises above **0.2 psi (A)** atmospheric.


1. If the verification is being carried out with another process (i.e. fitting the Magnet Service Syphon (except for Mobile systems), or replacing the Cold head), follow the **re-pressurization** part of this procedure when those processes are complete.

- **At Real-Time Values** - Monitor the **magnet pressure**.
  - Select **Update** to refresh data until pressure reaches  $\geq 15.6$  psi (A).



Data is only accurate from the last refresh of the page, or from when the page was opened last.

Fig. 50: Check magnet pressure level

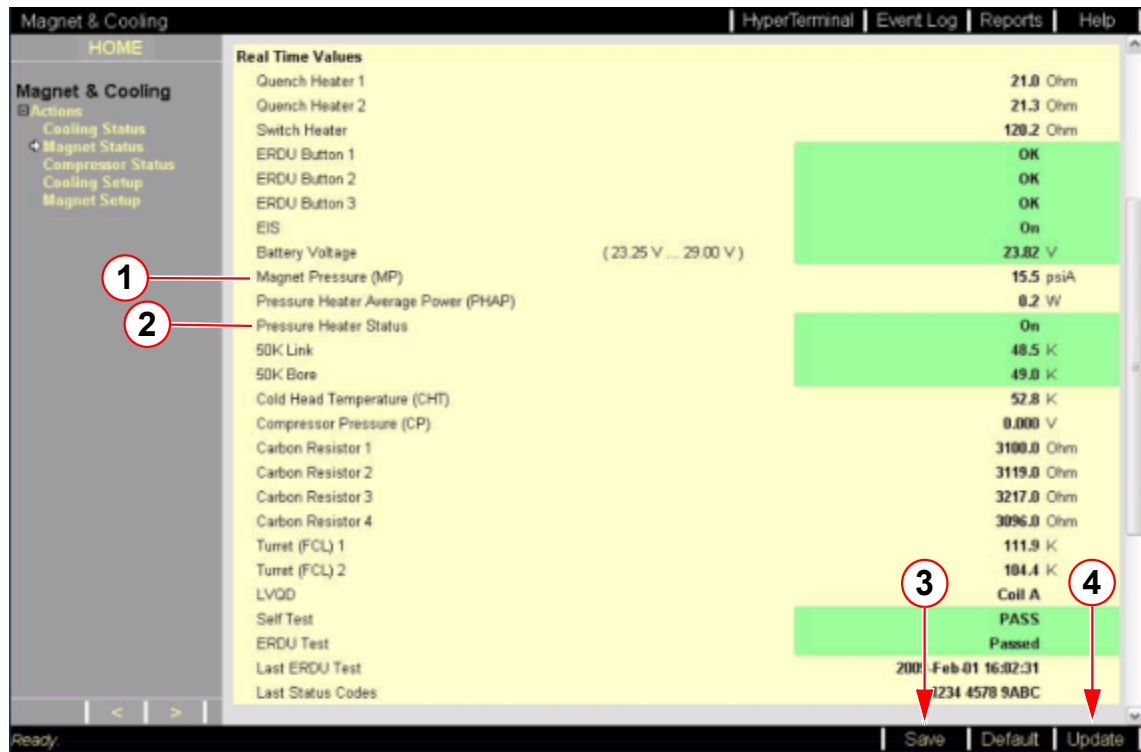

- (1) Magnet pressure (MP)
- (2) Pressure heater status
- (3) Save button
- (4) Update button

- **At MSUP Control** - Sets the pressure heater mode select the **Automatic mode** button (→ 1/ Fig. 47 Page 46)
- Click **UPDATE**. (→ 4/ Fig. 47 Page 46)
- Cancel any **alarms**.



Check that the bypass valve is locked securely in the closed position.



Contact the Magnet HSC if the test fails.



As the cold head recondenses the helium gas, the system pressure is returned to 15.5 psi (A).

## 2.3 Bottom filling the magnet



### CAUTION

Cold head failure.

The cold head seals can fail, and leak to atmosphere.

- » Always ensure the cold head bypass valve is **CLOSED** during helium filling.



The bottom fill procedure must **ONLY** be performed when the magnet is at **ZERO FIELD**.

**Bottom filling of the magnet will be required if:-**

- The liquid helium level -
  - Falls to **zero** percent,
  - **and**
  - The temperatures measured from the ceramic carbon resistors (**CCR**) and coil diodes, indicates temperatures **no warmer than** 110K. Refer to (→ Fig. 51 Page 51), (→ Fig. 52 Page 52), (→ Fig. 53 Page 53).

The magnet must be filled through the **magnet syphon cone** to allow filling from the bottom of the helium vessel.



Inform Magnet HSC immediately if the helium level is zero. Magnet coil temperatures must be confirmed before continuing with the helium fill.

**Possible causes for this include:-**

- Post magnet quench, (ramping or from persistence).
- Faulty helium level probes.
- Alarm box buttons still depressed.
- Pressure heater is left on.
- Long transportation and storage times.
- Cold head or compressor failure.

### 2.3.1 Checking magnet temperatures

Fig. 51: OR98 Magnet Acceptance Document (example)

| <b>SIEMENS</b>                                                                                                                          |                               | <b>Magnet Acceptance Documentation</b> |                      |
|-----------------------------------------------------------------------------------------------------------------------------------------|-------------------------------|----------------------------------------|----------------------|
|                                                                                                                                         |                               | <b>Installation Data</b>               |                      |
| <b>MR Magnet Technology Quality Regulations</b>                                                                                         |                               |                                        |                      |
| <b>Cryostat Number</b>                                                                                                                  | 21236477                      |                                        | <b>SMT Reference</b> |
|                                                                                                                                         |                               | 21245226 05.10                         |                      |
| <b>System Description</b>                                                                                                               |                               |                                        |                      |
| 1.495T Actively Shielded Superconducting Twin Turret Magnet System                                                                      |                               |                                        | OR98                 |
| Measured Magnet Centre (Offset)                                                                                                         | +1 (P) mm                     |                                        |                      |
| Gradient Coil Position                                                                                                                  | 39 mm                         |                                        |                      |
| Overall System Length                                                                                                                   | 1409 mm                       |                                        |                      |
| This system requires 12 Z2 rings and 2.38 Kg of iron to shim. This is a theoretical calculation that may change due to site conditions. |                               |                                        |                      |
| <b>Operating Parameters</b>                                                                                                             |                               |                                        |                      |
| <b>Description</b>                                                                                                                      | <b>Measured value</b>         |                                        |                      |
| Operating Current @ 1.495T                                                                                                              | 510.83 AMPS                   |                                        |                      |
| Field/Current Ratio                                                                                                                     | 12.46 kHz/0.1 Amp             |                                        |                      |
| Magnet Type: OR98-??                                                                                                                    | <b>Carbon Resistor Values</b> |                                        |                      |
| <b>POSITION</b>                                                                                                                         | <b>Room Temperature</b>       | <b>Helium Temperature</b>              |                      |
| CCR4 TDC Outer (D-E)                                                                                                                    | 926                           | 2951                                   |                      |
| CCR3 Top Inner (C-E)                                                                                                                    | 1060                          | 2812                                   |                      |
| CCR2 Middle Inner (B-E)                                                                                                                 | 1035                          | 2977                                   |                      |
| CCR1 Bottom Inner (A-E)                                                                                                                 | 1093                          | 3138                                   |                      |
| TDC Calibration Values:                                                                                                                 | 882.4                         | 1181.8                                 | 3075.3               |
| <b>Helium Level Probe Resistance (When using a 250 mA supply)</b>                                                                       |                               |                                        |                      |
|                                                                                                                                         | Probe 1                       | Probe 2                                |                      |
| 0% Resistance                                                                                                                           | 107.0                         | 107.0 Ohms                             |                      |
| 0 % Calibration                                                                                                                         | 448                           | 448                                    |                      |
| 100 % Resistance                                                                                                                        | 17.9                          | 17.9 Ohms                              |                      |
| 100 % Calibration                                                                                                                       | 75                            | 75                                     |                      |
| Diode Table No                                                                                                                          | -2                            |                                        |                      |
| Pressure Heater Res                                                                                                                     | 22.1                          | Ohms                                   |                      |
| <b>Final Temperature Sensor Recon Values</b>                                                                                            |                               |                                        |                      |
| 50Kbore                                                                                                                                 | 50K Link                      | Coldhead                               | Recon Margin         |
| 56                                                                                                                                      | 55                            | 41 K                                   | 782 Milliwatts       |
| Checked by                                                                                                                              | K Suraina                     |                                        | Signature            |
| Date                                                                                                                                    | 28/05/2010                    |                                        |                      |



Fig. 52: OR99 Magnet Acceptance Document (example)

|                                                                                                                                                       |                                                      |                                            |
|-------------------------------------------------------------------------------------------------------------------------------------------------------|------------------------------------------------------|--------------------------------------------|
| <b>SIEMENS</b>                                                                                                                                        | Magnet Acceptance Documentation<br>Installation Data |                                            |
| MR Magnet Technology Quality Regulations                                                                                                              |                                                      |                                            |
| <b>Cryostat Number</b>                                                                                                                                | <u>21240447</u>                                      | <b>SMT Reference</b> <u>21243544 05.10</u> |
| <b>System Description</b>                                                                                                                             |                                                      |                                            |
| <u>2.894T Actively Shielded Superconducting Twin Turret Magnet System</u>                                                                             |                                                      | <u>OR99</u>                                |
| Measured Magnet Centre (Offset)                                                                                                                       | <u>-3 (S)</u>                                        | mm                                         |
| Gradient Coil Position                                                                                                                                | <u>43</u>                                            | mm                                         |
| Overall System Length                                                                                                                                 | <u>1669</u>                                          | mm                                         |
| This system requires <u>12</u> Z2 rings and <u>1.44</u> Kg of iron to shim. This is a theoretical calculation that may change due to site conditions. |                                                      |                                            |
| <b>Operating Parameters</b>                                                                                                                           |                                                      |                                            |
| <u>Description</u>                                                                                                                                    | <u>Specification value</u>                           | <u>Measured value</u>                      |
| Operating Current @ 2.894T                                                                                                                            | < 500 AMPS                                           | <u>465.97</u> AMPS                         |
| Field/Current Ratio                                                                                                                                   | N/A                                                  | <u>26.44</u> kHz/0.1 Amp                   |
| Magnet Type: <u>OR99-??</u>                                                                                                                           |                                                      |                                            |
| <b>Carbon Resistor Values</b>                                                                                                                         |                                                      |                                            |
| <u>POSITION</u>                                                                                                                                       | <u>Room Temperature</u>                              | <u>Helium Temperature</u>                  |
| CCR4 TDC (D-E)                                                                                                                                        | <u>1072</u>                                          | <u>2863</u>                                |
| CCR3 Top (C-E)                                                                                                                                        | <u>872</u>                                           | <u>3114</u>                                |
| CCR2 Middle (B-E)                                                                                                                                     | <u>1027</u>                                          | <u>3182</u>                                |
| CCR1 Bottom (A-E)                                                                                                                                     | <u>1036</u>                                          | <u>2963</u>                                |
| TDC Calibration Values:                                                                                                                               | <u>1022.6</u> <u>1322.1</u>                          | <u>2991.4</u>                              |
| <b>Helium Level Probe Resistance (When using a 250 mA supply)</b>                                                                                     |                                                      |                                            |
|                                                                                                                                                       | <u>Probe 1</u>                                       | <u>Probe 2</u>                             |
| 0% Resistance                                                                                                                                         | <u>80.3</u>                                          | <u>80.3</u> Ohms                           |
| 0 % Calibration                                                                                                                                       | <u>336</u>                                           | <u>336</u>                                 |
| 100 % Resistance                                                                                                                                      | <u>0.5</u>                                           | <u>0.5</u> Ohms                            |
| 100 % Calibration                                                                                                                                     | <u>002</u>                                           | <u>002</u>                                 |
| Diode Table No                                                                                                                                        | <u>-2</u>                                            |                                            |
| Pressure Heater Res                                                                                                                                   | <u>20.5</u> Ohms                                     |                                            |
| <b>Final Temperature Sensor Recon Values</b>                                                                                                          |                                                      |                                            |
| <u>50Kbore</u>                                                                                                                                        | <u>50K Link</u>                                      | <u>Coldhead</u>                            |
| <u>54</u>                                                                                                                                             | <u>50</u>                                            | <u>40</u> K                                |
|                                                                                                                                                       |                                                      | <u>Recon Margin</u>                        |
|                                                                                                                                                       |                                                      | <u>874</u> Milliwatts                      |
| Checked by                                                                                                                                            | <u>K Suraina</u>                                     | Signature <u></u>                          |
| Date                                                                                                                                                  | <u>10/05/2010</u>                                    |                                            |
| Page 1 of 1                                                                                                                                           |                                                      | 25-AD EPS 000 01                           |



Fig. 53: OR103 Magnet Acceptance Document

| SIEMENS                                                                    |                                                                          | Magnet Acceptance Documentation          |                           |
|----------------------------------------------------------------------------|--------------------------------------------------------------------------|------------------------------------------|---------------------------|
|                                                                            |                                                                          | Installation data                        |                           |
|                                                                            |                                                                          | MR Magnet Technology Quality Regulations |                           |
| <b>Cryostat number</b>                                                     | 33002284                                                                 | <b>SMT Shipped System</b>                | 33002281 07.11            |
| <b>System description</b>                                                  | 2.894T Actively Shielded Superconducting Twin Turret Magnet System OR103 |                                          |                           |
|                                                                            | <b>SKYRA</b>                                                             | <b>MONET</b>                             |                           |
| <b>Measured magnet centre (Offset)</b>                                     | -4 (S)                                                                   | -4 (S)                                   | mm                        |
| <b>Gradient Coil Position</b>                                              | 174                                                                      | 174                                      | mm                        |
| <b>Overall Magnet Length</b>                                               | 1670                                                                     | 1670                                     | mm                        |
| <b>This system requires</b>                                                | 0                                                                        | 0                                        | Z2 Rings                  |
| <b>This system requires</b>                                                | 2.63                                                                     | 2.63                                     | kg of iron to shim        |
| (This is a theoretical calculation that may change due to site conditions) |                                                                          |                                          |                           |
| <b>Operating parameters</b>                                                |                                                                          |                                          |                           |
| <b>Description</b>                                                         | <b>Specification value</b>                                               | <b>Measured value</b>                    |                           |
| Operating current @ 2.894T                                                 | < 500 Amps                                                               | 487.35 Amps                              |                           |
| Field/Current Ratio                                                        | N/A                                                                      | 25.28 kHz/0.1 Amp                        |                           |
| Magnet type: OR103-00                                                      |                                                                          |                                          |                           |
| <b>Carbon resistor values</b>                                              |                                                                          |                                          |                           |
| <b>Position</b>                                                            | <b>Room temperature</b>                                                  | <b>Nitrogen temperature</b>              | <b>Helium temperature</b> |
| CCR4 d - e top                                                             | 1071 Ohms                                                                |                                          | 2573 Ohms                 |
| CCR3 c - e                                                                 | 991 Ohms                                                                 |                                          | 2612 Ohms                 |
| CCR2 b - e                                                                 | 929 Ohms                                                                 |                                          | 2797 Ohms                 |
| CCR1 a - e bottom                                                          | 1007 Ohms                                                                |                                          | 2841 Ohms                 |
| CCR4 calibration values                                                    | 1031 Ohms                                                                | 1304 Ohms                                | 2621 Ohms                 |
| <b>Helium level probe resistance (when using a 250mA supply)</b>           |                                                                          |                                          |                           |
|                                                                            | <b>Probe 1</b>                                                           | <b>Probe 2</b>                           |                           |
| 0% resistance                                                              | 107 Ohms                                                                 | 107 Ohms                                 |                           |
| 0% calibration                                                             | 448                                                                      | 448                                      |                           |
| 100% resistance                                                            | 10.8 Ohms                                                                | 10.8 Ohms                                |                           |
| 100% calibration                                                           | 45                                                                       | 45                                       |                           |
| Diode table no.                                                            | 2                                                                        |                                          |                           |
| Pressure heater resistance                                                 | 20.5 Ohms                                                                |                                          |                           |
| <b>Final temperature sensor recon values</b>                               |                                                                          |                                          |                           |
| Shield bore                                                                | Shield Link                                                              | Coldhead                                 | Recon margin              |
| 55K                                                                        | 48K                                                                      | 39K                                      | 853mW                     |
| <b>Checked by</b> K. Baker                                                 | <b>Signed</b> <i>Km. Baker</i>                                           |                                          |                           |
| <b>Date</b> 11/07/2011                                                     |                                                                          |                                          |                           |
| Page 1 of 1                                                                |                                                                          | 25-AD EPS 000 01                         |                           |

To **recover** a magnet with a **bottom helium fill**, the following must be checked:-

Access the Service software at the **MRC (MRAWP) Host** to measure:

- **Helium level** values.
- **CCR** values.

Tab. 5 Expected CCR resistor values

| 41 way plug, measured in Ohms       | Expected values room temperature | Expected values 77 K temperature | Expected values 4.2 K temperature |
|-------------------------------------|----------------------------------|----------------------------------|-----------------------------------|
| All ceramic carbon resistors (CCRs) | 880 - 1050                       | 1150 - 1350                      | 2800 - 3380                       |



Carbon 4/CCR4 is for measuring helium gas temperatures during ramping and is not applicable to filling.

### 2.3.2 Bottom fill procedure




The bottom fill procedure must **ONLY** be performed when the magnet is at ZERO FIELD.

- Prepare the **dewar** and **dewar syphon**. (→ Dewar preparation / Page 15), (→ Dewar syphon preparation / Page 17)
- Check the **ceiling height**. (→ Low ceiling height / Page 41)
  - Assemble the **syphon** (only if it can be fitted), or
  - Assemble the **diffuser**, bottom **extension** and support **clamp**.

### 2.3.2.1 Assemble the service syphon

Fig. 54: Magnet service syphon assembly - bottom filling

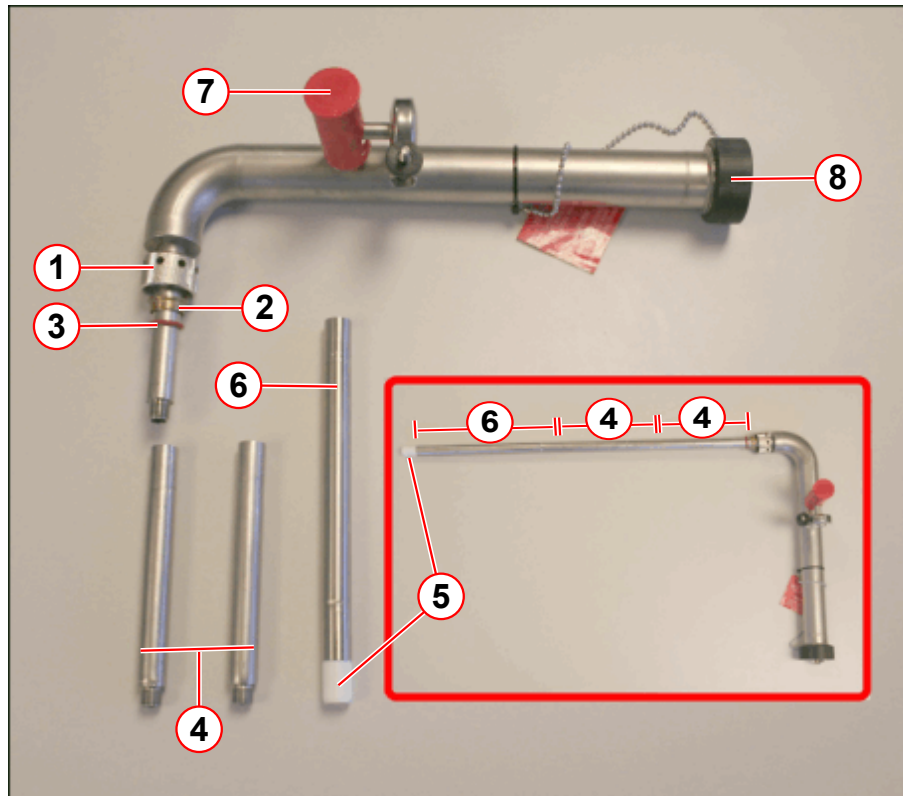


- (1) Ring nut
- (2) Brass spacer
- (3) O-ring seal
- (4) 140 mm syphon legs (2 off)
- (5) 22 mm PTFE seal (bottom fill)
- (6) 200 mm syphon leg (1 off)
- (7) Vacuum valve
- (8) Syphon end cap

Tab. 6 OR98 / OR99 / OR103 BOTTOM filling syphon extension legs

| Magnet type     | 140 mm leg | 200 mm leg | 22 mm PTFE seal |
|-----------------|------------|------------|-----------------|
| OR98 (Aera)     | 2          | 1          | 1               |
| OR99/103 (Skra) |            |            |                 |

Fig. 55: Assembling the Ring Nut Seal

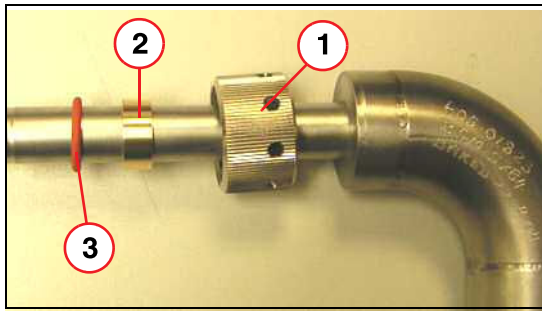

- (1) Ring nut
- (2) Brass spacer
- (3) Silicon O-ring

- Fit the ring nut, brass spacer and the O-ring to the **service syphon leg**.

Refer to table above (→ Tab. 6 Page 55)

- Connect the **extension legs** and the 22 mm **PTFE seal**:
  - Screw the 2 x **140 mm legs** together, and connect to the **syphon**.
  - Fit the 1 x **200 mm leg**.
  - Fit the **22 mm PTFE seal** to the 200mm leg.



The minimum length required to reach the syphon cone is:


OR98 - 554 mm (+/- 5 mm)

OR99 - 540 mm (+/- 5 mm)

### 2.3.2.2 Depressurize the helium vessel

- **De-activate** the **pressure heater**. (→ Deactivate the pressure heater / Page 21)
- **Depressurize** the helium vessel (→ Depressurize the helium vessel / Page 22)
- Note the **pressure** (using the **MRC**, or **Laptop**) when the **compressor** switches **OFF**.

Fig. 56: WIKKA Pressure switch - with a small positive pressure

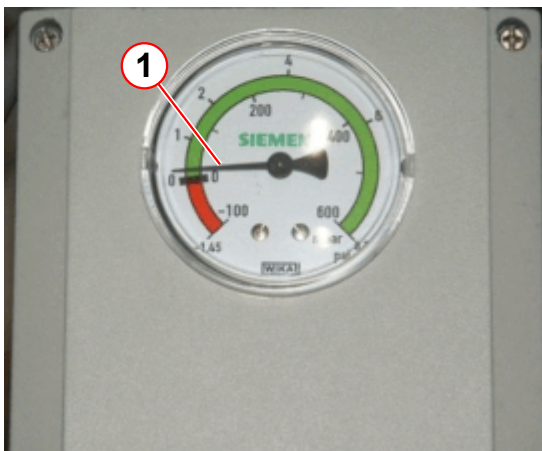

- **Check** the helium vessel has a **small positive pressure**.

- **Warm** and dry the **syphon blanking plug** using a heat gun to ease removal.

- Dry the syphon port.
- Remove the syphon **end cap**. (→ Fig. 33 Page 33)
- Remove the syphon **blanking plug**.



### CAUTION


Risk of cold burns.

Cold gas will exhaust from the syphon port on opening.

- » Stand clear when opening.
- » Wear protective gloves and safety glasses.


- Slowly lower the **service syphon** fully into the neck.
- The **syphon leg** must locate fully into the **syphon cone** - make sure it is **fully engaged**.

Fig. 57: Bottom fill syphon cone (pictorial view)

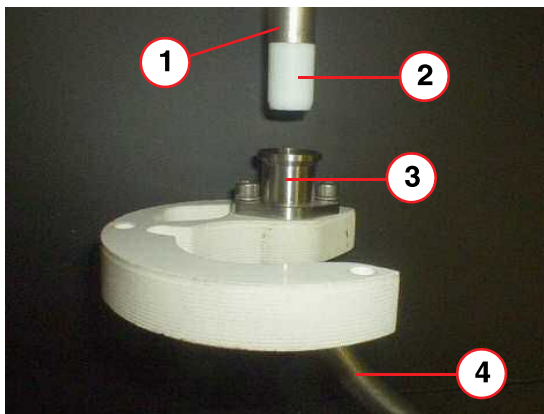

- (1) Syphon leg
- (2) PTFE Seal (Bottom Fill)
- (3) Syphon cone
- (4) Syphon tube

Fig. 58: Bottom fill - syphon leg fully engaged (pictorial view)

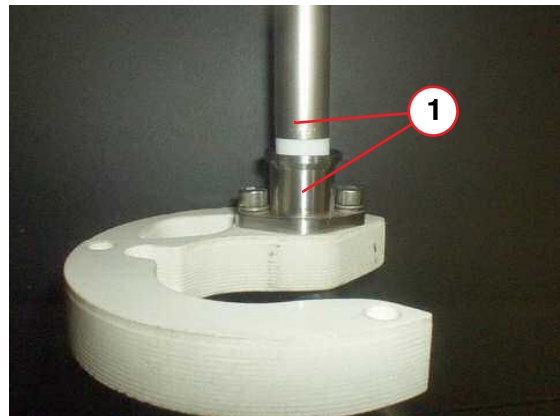

- (1) Syphon leg - fully engaged



Helium efficiency will be reduced if the syphon leg is not fully engaged in the syphon cone.

- Tighten the **syphon port ring nut** (→ 1/Fig. 55 Page 56) by hand, then apply a half turn, (using the appropriate spanner), to form an air-tight seal.
- Re-fit the **syphon end cap**.

### 2.3.3 Liquid helium transfer - bottom fill



Before helium transfer, refer to (→ Helium pressures at altitude / Page 42)




The helium dewar must be positioned close to the magnet, with the dewar syphon fitted.



For bottom filling a magnet at zero, there is no need to pre-cool both syphons.

- Check the dewar gas supply pressure is set at **0.27bar (4 psi (G))**.
- Open the **dewar syphon flow valve** half way (2 turns) to **purge** the syphon of air.

Fig. 59: Dewar syphon flow valves - open and closed positions

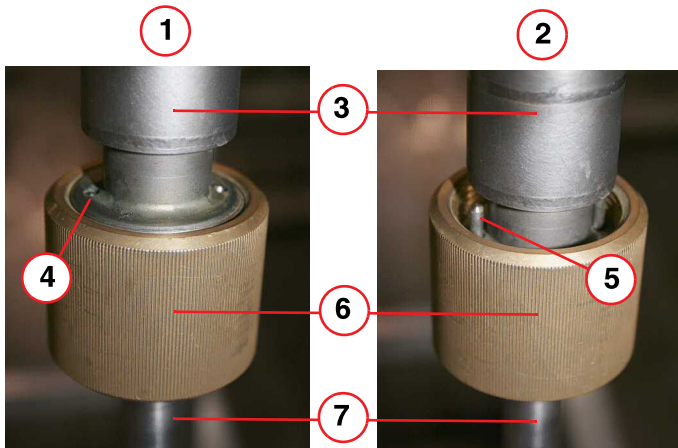

- (1) FULLY OPEN position
- (2) FULLY CLOSED position
- (3) Stainless steel body
- (4) Guide pin flush with body
- (5) Guide pin
- (6) Brass nut
- (7) Syphon leg

Fig. 60: Helium filling



- Close the dewar syphon **flow valve**.
- Remove the service syphon end cap and **couple the two syphons together** and secure the **locknut**.
- Open the dewar syphon **flow valve**.
- Maximize the transfer by fully opening the **dewar syphon flow valve**.
  - » The flow valve is **fully open** when the two **guide pins** at the top of the brass nut are flush with the top of the stainless steel body. (→ 4/ Fig. 59 Page 58)
  - » From **fully closed** to **fully open**, is **4 turns**.
- Ensure the transfer dewar pressure is **4 psi (A)** at normal altitude. Refer to (→ Helium pressures at altitude / Page 42).

### 2.3.3.1 Monitor the CCRs

- Monitor the **CCR**s every **15 minutes**, to measure the **temperature drop**.
  - » Refer to the **acceptance documentation** (Installation Protocol) for required values.
  - » If the CCRs **do not** get cooler, **contact Magnet HSC**.
- When the required value is reached, **monitor** the **helium level** at the **Host** or **MSUP**.

### 2.3.3.2 Perform a helium level check



During installation, if the MRC is not available:

1. Take a helium sample using the MSUP. (→ Helium level check at the MSUP / Page 61)


#### Helium level check at the MRC

- Check the liquid helium level by using Service Software on MRC (MRAWP) Host.
  - Check the liquid helium level using the MRC.
  - Access the Service Software at the MRC.
  - Select the menu item Magnet & Cooling - Magnet Status.



- Press the Start button (→ 3/ Fig. 61 Page 60) to initiate the helium measurement.  
(→ 1/ Fig. 61 Page 60)

Fig. 61: Checking liquid helium and magnet pressure levels

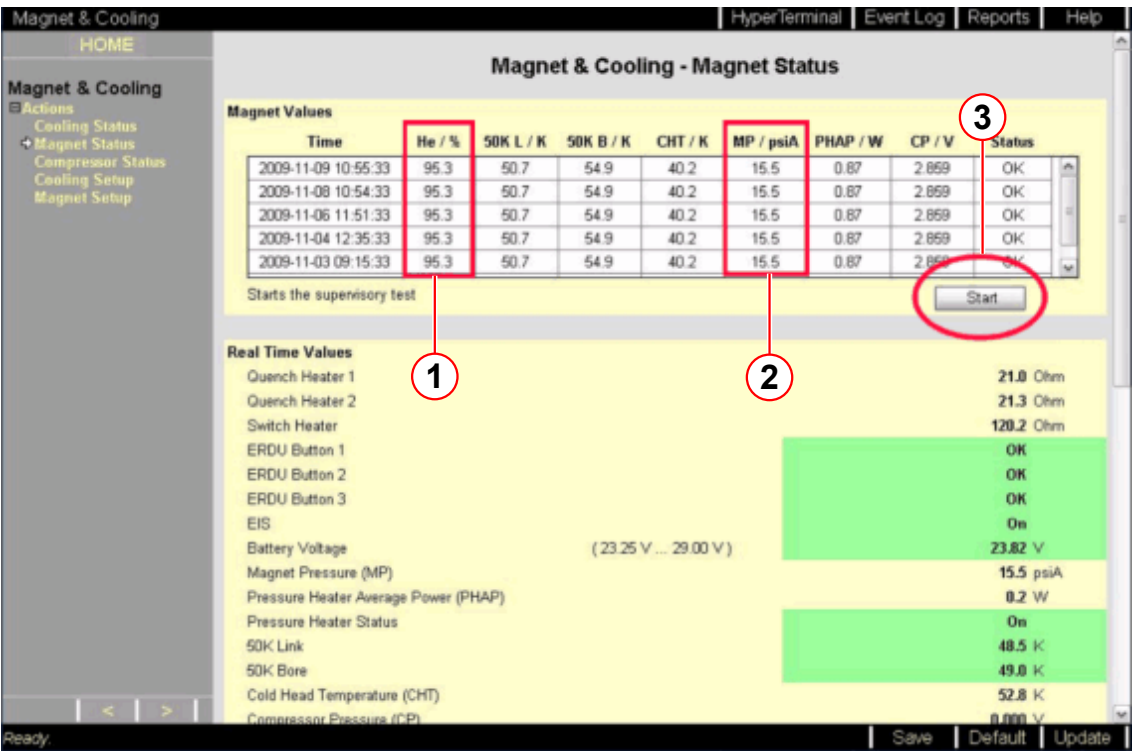


- (1) Liquid helium level
- (2) Magnet pressure
- (3) Start button

**i** Data is only accurate from the last refresh of the page, or from when the page was opened last.


## Helium level check at the MSUP

Fig. 62: MSUP 2306 - depress SAMPLE to activate

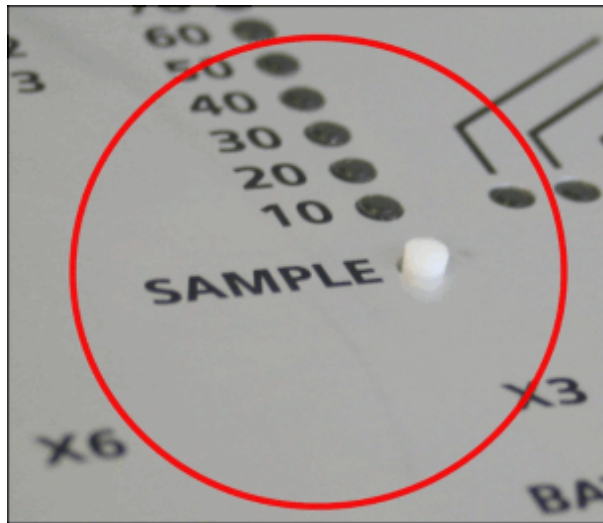

- Locate the SAMPLE switch on the MSUP lid.
- Depress the sample switch for 2 s, and release.

Fig. 63: MSUP - Helium level sample

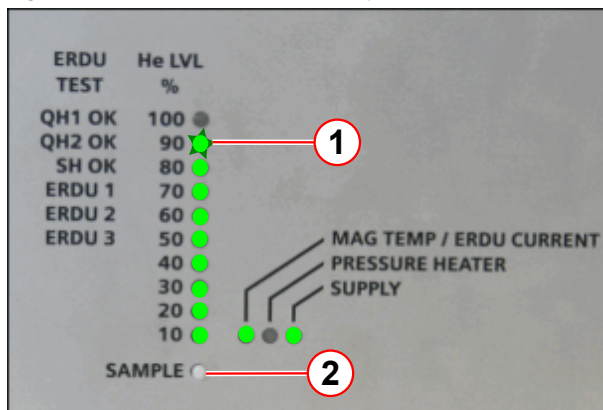

- (1) Sample switch  
(2) Helium level

- The green LEDs illuminate and display the actual liquid helium level percentage.
- If the 90% LED is flashing, STOP the transfer. Refer to (→ Tab. 7 Page 61).

Tab. 7 MSUP Helium display

| MSUP LEDs                      | Helium level | Comments or actions            |
|--------------------------------|--------------|--------------------------------|
| 0 - 80% ON                     | As displayed | Continue to fill as required.  |
| 90% ON                         | As displayed | Stop transfer.                 |
| *90% FLASHING<br>10 - 80% ON   | 95% - 97%    | STOP the transfer IMMEDIATELY! |
| *100% FLASHING<br>10% - 90% ON | 97% and over | Magnet is OVERFILLED.          |



\* The LEDs return to ON after approximately 30 seconds.



To avoid overfilling, the helium transfer must be stopped as soon as the 90% LED illuminates.

- Repeat the measurement every 5 minutes. An increase in level should be seen within 20 minutes.
  - » If not, stop the transfer and investigate the reason for poor transfer efficiency.



Opening the bypass valve will help to reduce the back pressure, enabling more efficient helium transfer.



Maintain a pressure of about 0.27 bar (4 psi (G)) at the helium dewar. Too much or too little pressure will result in inefficiency.

At 4 psi (G), the transfer rate is approximately 10 liters per minute. **Do not leave the transfer unattended.** A 250 liter dewar will empty in 25 - 30 minutes.

If frost or condensation forms on the body of either syphon, it is an indication of poor vacuum. **Stop the transfer immediately - the syphon may need to be replaced.**

## NOTICE

**Risk of quench.**

**Warm helium gas will cause a quench.**

- » **Stop the transfer immediately in case you hear whistling (a high pitched sound) from the syphon, or the dewar pressure falls.**

### 2.3.4 Stopping the helium transfer

Fig. 64: Dewar manometer



- As soon as the required helium level is obtained or a maximum of level of **95%**, or the dewar is empty, stop the transfer by closing the **flow valve** at the dewar syphon. **Do not fill above 95%**. (→ Helium level check at the MRC / Page 59)
  - Whistling (high pitched sound) of the syphon or noticeable drop in dewar pressure are indications of an empty dewar.
- The **dewar syphon flow valve** must be closed **immediately** the dewar empties. (→ 6/ Fig. 59 Page 58)
- Set the **dewar** gas pressure regulator to **zero**.
- Close the dewar **primary ball valve**. (→ 1/ Fig. 14 Page 19)
- **Remove** the **gas line** from the **dewar**.
- Refer to (→ Changing supply dewars during transfer / Page 36) if the magnet is not full.

#### NOTICE

Risk of quench.

Warm helium gas will cause a quench.

- » Stop the transfer immediately in case you hear whistling (a high pitched sound) from the syphon, or the dewar pressure falls.

### 2.3.5 Removing the service syphon

- **Deactivate** the pressure heater. (→ Deactivate the pressure heater / Page 21)
- **Depressurize** the helium vessel. (→ Depressurize the helium vessel / Page 22)

### ⚠ CAUTION

Danger of air ingress.

Internal ice blockage could cause excessive pressure.

- » Warm up syphon port plug sufficiently to ensure there is no ice on the threads when tightening.

- Warm the **syphon joint** to room temperature (using a heat gun).
- Warm the syphon **entry valve** to room temperature and **dry**.
- Warm the syphon **end cap** ready for fitment onto the syphon.
- Disconnect the syphons.

### ⚠ CAUTION


Risk of cold burns.

Cold gas will exhaust from the syphon port on opening.

- » Stand clear when opening.
- » Wear protective gloves and safety glasses.


- Immediately close off the **end** of the service syphon with the **syphon end cap**.
- Immediately **remove** the service **syphon** from the magnet.

Fig. 65: Syphon blanking plug

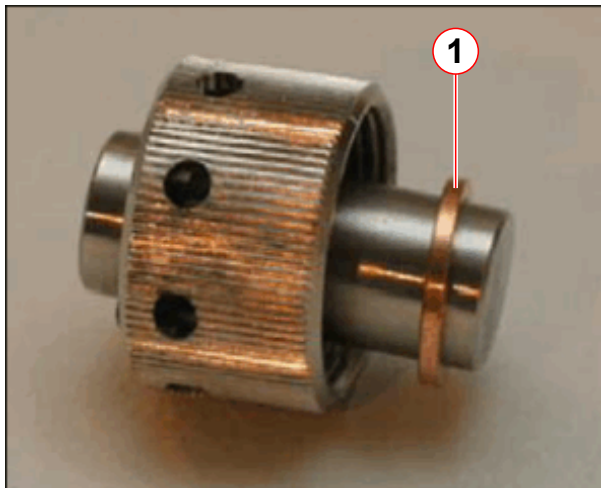

(1) Copper seal

Fig. 66: Fit/remove the syphon blanking plug

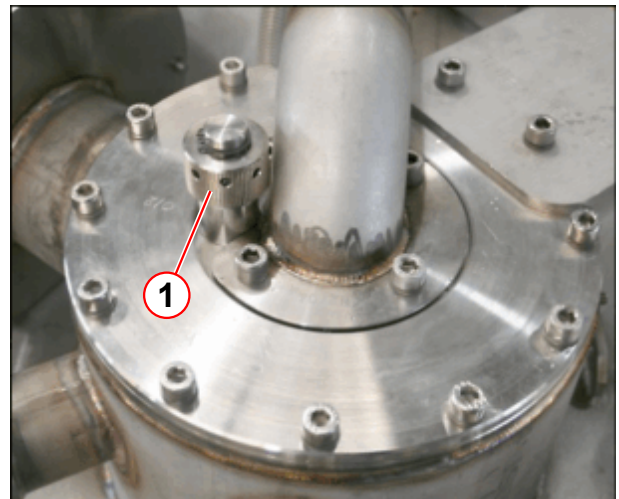

- Immediately fit the **syphon blanking plug** using a new copper seal (→ Fig. 65 Page 64).
- Warm up the threads on the plug and **carefully** tighten with the syphon C-spanner.
  - There might be **ice** on the threads - **do not use undue force**. Damaged threads will result in a leak.
- Warm the **top plate** and **syphon plug** using a heat gun. **Re-tighten** the plug.

- Warm up the magnet bypass **valve assembly** and check the **bypass valve** is closed.

**CAUTION**

Danger of air ingress.

Internal ice blockage could cause excessive pressure.

- » Warm up syphon port plug sufficiently to ensure there is no ice on the threads when tightening.

- Pressurize the helium vessel to 15.5 psi (A).
  - Reactivate the pressure heater.(→ Reactivate the pressure heater / Page 47)
  - Check that all **pressure alarms** are deactivated.
  - Check the **compressor** is running.
- Warm the **extension legs** using a heat gun.
- Remove the **bottom fill extension leg** and **PTFE seal**.
- Warm the **turret and venting valves** to room temperature with a heat gun.
  - remove any **surplus water** with paper towelling.
- Move the dewar out of the imaging suite into a **well ventilated area**.
- Complete **depressurizing** the dewar by opening the **primary ball valve**. (→ 5/Fig. 38 Page 37)
- When at the gas pressure is at **zero**, close the **primary ball valve** on the dewar.
- Remove the dewar syphon.
- Immediately close the **syphon entry valve** on top of the dewar. (→ 4/Fig. 38 Page 37)
- Open the **low-pressure dewar relief valve**. (→ 1/Fig. 38 Page 37)
- Warm both **syphons** with the heat gun **until dry**, and store safely.
- Leak check the turret and venting.




After a bottom fill, it is recommended that the first magnet energization is performed the next day. This will allow the magnet coils to thermally stabilize.

### 2.3.6 Final work steps

After transfer, check:-


- Press the **LVQD reset button** to ensure the functionality of the LVQD circuit.

Fig. 67: LVQD reset button

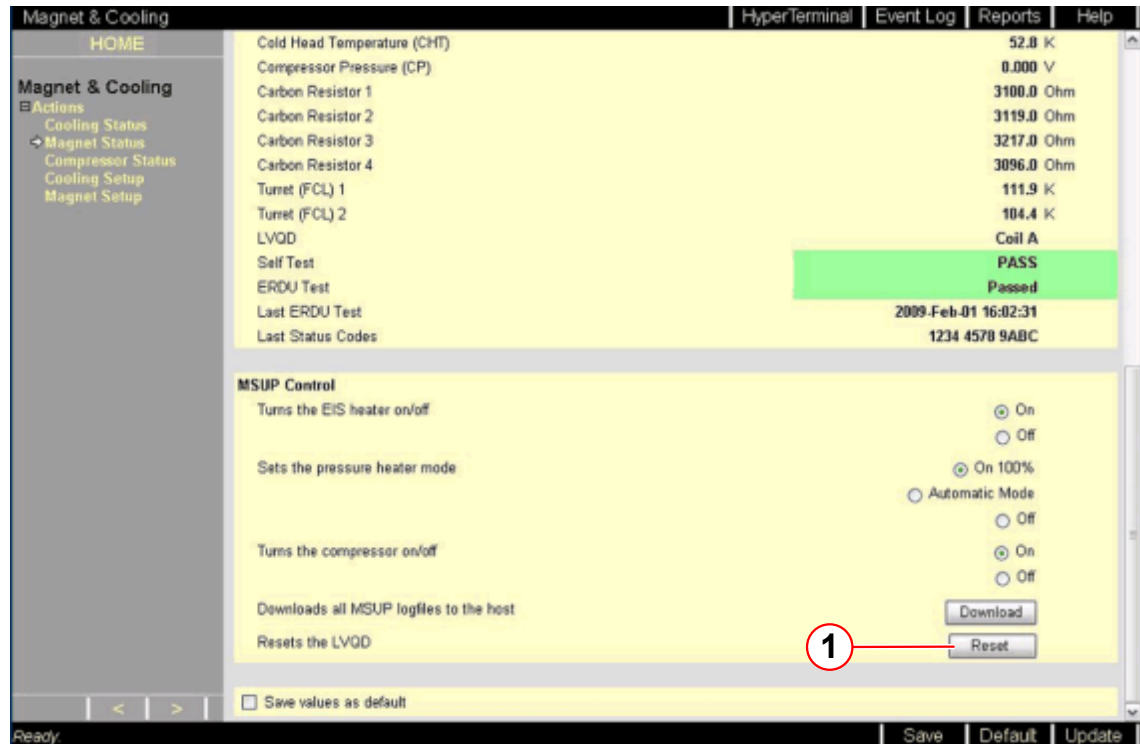

- The magnet bypass valve is closed.
- The pressure heater is operating correctly.

Fig. 68: MSUP - Pressure heater LED

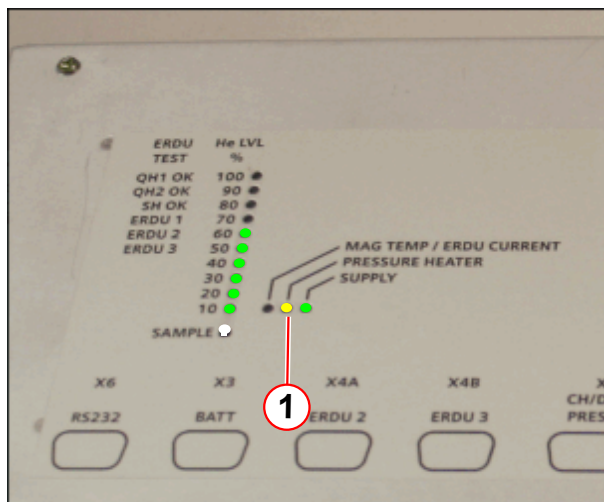

(1) Pressure heater LED - yellow

- No alarm signals are evident at the MSUP or the alarm box.
- The compressor is switched on and running.
- The magnet pressure is stable at 15.5psi (A).

- The pressure heater LED on the MSUP is illuminated, or switching on and off.



- The **turret** and **venting** system are **warm** and **dry**.
- The filling **port** and **venting** system have **no leaks**. (Refer to **Venting Leak Check procedure**).
- The 4K **cold head flange** has **no leaks**.
- The used **dewars** are correctly **closed**, with their low pressure relief valves **open**.

### 3.1 Quench tube check

- Perform a visual inspection of the quench tube.
  - » Refer to (→ Quench tube / M7-010.833)
- Check if the Installation protocol is signed by the Quench line manufacturer.



Ensure the Quench Line installation is complete before performing the first ramping procedure.

The Quench Line Acceptance Protocol must be checked before the first ramp.

Check that the Quench Line Acceptance Protocol is signed off by the Quench Line Manufacturer.

Acceptance Protocol can be found under:

**Product Information Home > Planning > Planning Guides > MR Systems > Quench Line Design and Calculation Tool**

### 3.2 Identification of the controlled access area (0.5 mT area)

- Check that the warning signs identifying the controlled access area (0.5 mT zone) have been posted at the entrance to the RF room as well as all other areas leading to the 0.5 mT zone.
  - » Refer to: (→ Identification of the controlled access area (0.5 mT area) / M7-010.833.05)



#### **WARNING**

##### **Strong magnetic field**

Failure to observe the following safety warning may result in physical injury and equipment damage.

- » Before ramping up the magnet, ensure that no ferrous tooling or equipment are present in the imaging suite.
- » Ensure that the warning signs identifying the controlled access area (0.5 mT zone) have been posted at the entrance to the RF room as well as all other areas leading to the 0.5 mT zone.
- » The sign must be posted in a way that it is visible even when the door to the controlled area is opened.

### 3.3 Magnet stop button test

- Perform the Magnet stop button test
  - » Refer to (→ ERDU / M7-010.833)

## 3.4 Venting leak check

### 3.4.1 General

For all recondensing systems, it is essential that the magnet system is leak tight to maintain a zero-loss cryogenic performance.

Individual leaks must be  $1 \times 10^{-4}$  mbar l/s or less.

Leak checking must be performed with an electronic leak detector capable of detecting leaks, as small as  $1 \times 10^{-4}$  mbar.

The magnet must be leak checked after a service operation has been performed, such as:

- Installation and Startup.
- After helium filling.
- After replacing a component on the high pressure side of the magnet. (cold head, venting).
- After a magnet quench. (Post quench).



#### CAUTION

##### Risk of icing


Small leaks may not show any change in pressure. Over time, they will contribute to air ingress and icing.

- » Leak check with the recommended equipment for small leaks greater than  $1 \times 10^{-4}$ .

### 3.4.2 Venting System

The document describes the procedure for **OR98**, **OR99** and **OR103** venting systems.

Fig. 69: Venting configuration

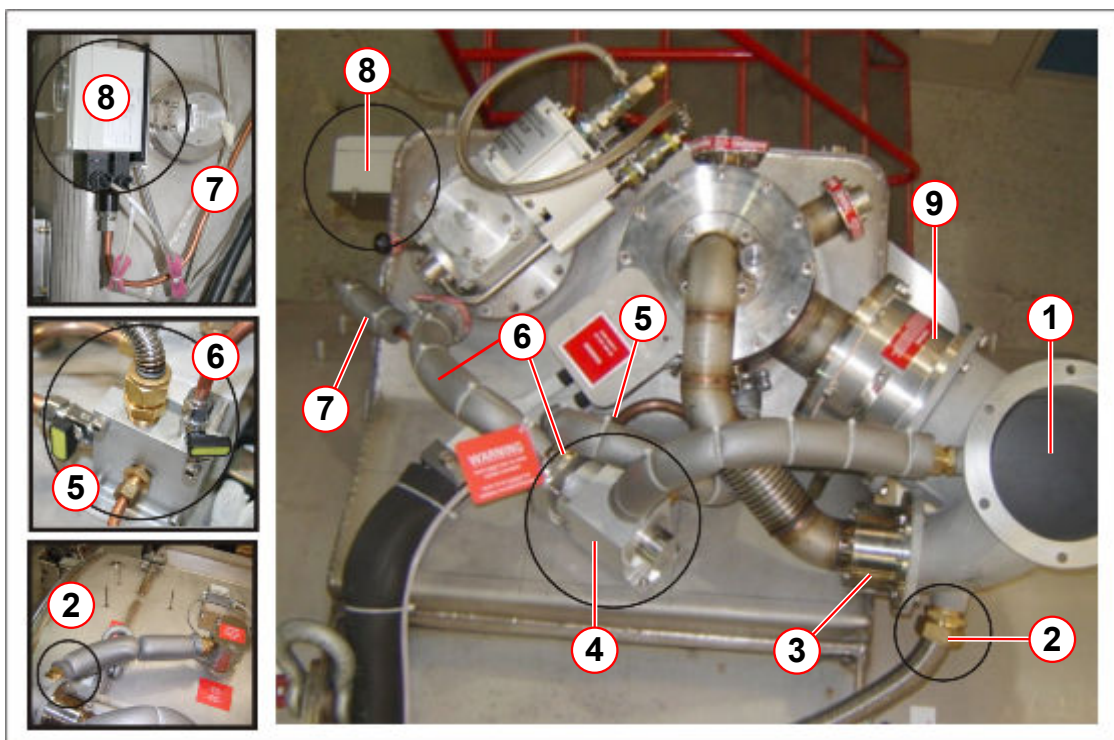

- (1) Quench valve elbow
- (2) OVC vent burst disk (if fitted)
- (3) Auxiliary burst disk (1.38 bar)
- (4) Venting valve assembly
- (5) Magnet system bypass line
- (6) Cold head vent line and bypass valve
- (7) Pressure display & control switch line
- (8) Pressure display & control switch
- (9) Quench valve with integral 15 psi (G) burst disk



The OVC burst disk and DOP port vent line (Pos.2) are no longer fitted.

### 3.4.3 Safety



#### CAUTION

Pressurized helium vessel.

Cold helium gas.

Failure to observe the following safety information may result in physical injury.

- » Fittings on the high-pressure side must not be disconnected without first depressurizing the helium vessel.

**WARNING**

Magnet at field.

Personal injury or damage to equipment may occur.

- » Do not use magnetic tools for repairing or tightening seals when the magnet is at field.
- » Only non-magnetic equipment may be taken into the magnet room.

### 3.4.4 Tooling and Spare Parts

#### 3.4.4.1 Tooling

- Heat gun.
- Paper towelling.
- Standard Service hand tools.
- Leak detection equipment.

Fig. 70: Edwards Leak detector (example)



- Use the Edwards hand-held leak detector (material number **10614850**). This is the recommended equipment for leaks greater than  $1 \times 10^{-4}$

**NOTE:** The hand-held device is calibrated in ml/sec:-

$1 \times 10^{-4}$  mbar l/s =  $1 \times 10^{-4}$  ml/sec at normal pressure and temperature.

Refer to the equipment manufacturer's operating instructions.



The function of the Edwards leak detector can be compromised by magnetic fields if placed near the magnet bore.

#### 3.4.4.2 Spare Parts

Helium leaks from the magnet pressure vessel are repaired by replacing a damaged seal or failed component.

- A number of seals are supplied in the spares kit. These must be replaced when used on the magnet.



- Venting components will require ordering from CSML.

### 3.4.5 Leak Check Process at Installation and Service

#### 3.4.5.1 Installation

- Check the system pressure on arrival. If the system has low or zero pressure then a gross leak has occurred during shipment.
- Zero pressure indicates a gross leak. Immediately check the burst disks.
- Check any joints (clamps and seals) that have been broken or damaged during transportation.
- If magnet pressure is still low, check for leaks through the magnet 17 psi (A) vent valve.

#### 3.4.5.2 Service

- The turret, vent valve assembly and vent lines must be warmed up to room temperature after any service procedures have been performed, prior to leak check.
- Leak check all seals that have been subjected to low temperatures and the subsequent build up of frost. For example, magnets that have been helium filled, depressurized, or the magnet ramped.
  - Check joints broken and re-made during servicing. Replace damaged seals.
  - Check that all clamps and seals are tight. Do not over-tighten them, it may cause a leak, or make an existing leak worse.

#### 3.4.5.3 Pre-requisites

- The helium vessel must be at normal operating pressure (15.5 psi (A)) prior to leak checking.
- Activate the pressure heater (if required) to raise pressure.
- If pressure heater activation continues for more than 30 min with NO increase in the magnet pressure:-
  - Stop heater activation and check for any obvious gross leaks such as the syphon ring nut or burst disks.
- All joints must be dry and at room temperature before starting a leak check.
  - » Warm all components up using a heat gun and dry. Moisture can cause false readings.

#### 3.4.5.4 Leak Check Process

- Switch on the leak detector and allow it to warm up (refer to manufacturer's instructions).
- Check the equipment is working correctly (detecting helium). Check equipment before and after use.
  - Set the leak detector to its maximum sensitivity.

- Check the response time of the detector.
- The background helium level in the magnet room must not be saturated with helium gas.
  - » **DO NOT** leak check immediately after a quench or helium filling.



The background must be at least one order of magnitude lower than the smallest leak expected, (i.e. leaks of  $10^{-4}$  mbar l/s, require background of order  $10^{-5}$  mbar l/s).

- Start at the lowest point.
- Hold the tip of the leak detector 5 mm above the surface of the component. **Keep the tip above water condensation.**
- Slowly move the detector around the seals or clamped joints.
- Hold the detector at a leak point until the reading has stabilized.
- Move the detector away and check that the indication disappears.
- Repeat the process to check that a leak is confirmed.
- Complete leak check of turret and venting.
- Fix all leaks (TSG) and repeat leak check of the component and seal.



Do not tighten any seals when the components are cold and covered in frost. Always warm up with a heat gun and dry with paper towelling.

#### 3.4.5.5 Leak Check Requirements

The following points must be observed at all times:

- Before dismantling a joint, check that a replacement seal is available.
- Always **CLOSE** the magnet bypass valve as soon as the system is depressurized.
- Always plug open venting with paper towelling to prevent air ingress.
- **Do not** leave 2 joints open at the same time.
- **Do not** leave the system open for longer than is necessary to complete the repair.
- Check that the sealing surfaces are clean, dry, and undamaged.
- **Do not** replace metal seals with rubber or silicon seals.



**Do not** use any form of vacuum grease on O-rings and seals. Vacuum grease dries over time and may cause a seal to leak.

For all leaks, first try tightening the clamp or the securing screws.

- **Do not** overtighten clamps or screws.
  - » If tightening does not solve the leak, disassemble the joint and check for any damage to the sealing surfaces.
  - » Replace the high load seal.

- **DO NOT** re-use metal seals. Assemble the joint using a new seal. Ensure that the seal is undamaged and clean, and is seated correctly.
- After repair, leak test the entire venting system.

### 3.4.6 High Pressure Turret and Venting Components

The venting and components on the high pressure side of the magnet turret require leak checking to ensure that the re-condensing performance is achieved. This includes all connections too:

- Cold Head
- Turret top plate, including syphon port and venting.
- 16 or 17 psi (A) Vent valve assembly.
- Pressure gauge assembly.
- Vent line connections.



If no external leak is found but the system pressure is low, then it is possible that one of the internal valves is leaking. This can be caused by a cracked or broken burst disk, or a broken O-ring.




Before a leak check, the venting components must be warmed to room temperature, and dried with paper towelling.

#### 3.4.6.1 Cold Head

Leak check the cold head components:

- Check the hand-valve and cold head vent line seals at both ends.
- Cold head O-ring seals.
- Temperature sensor connector seals.


Fig. 71: Cold head valve closed

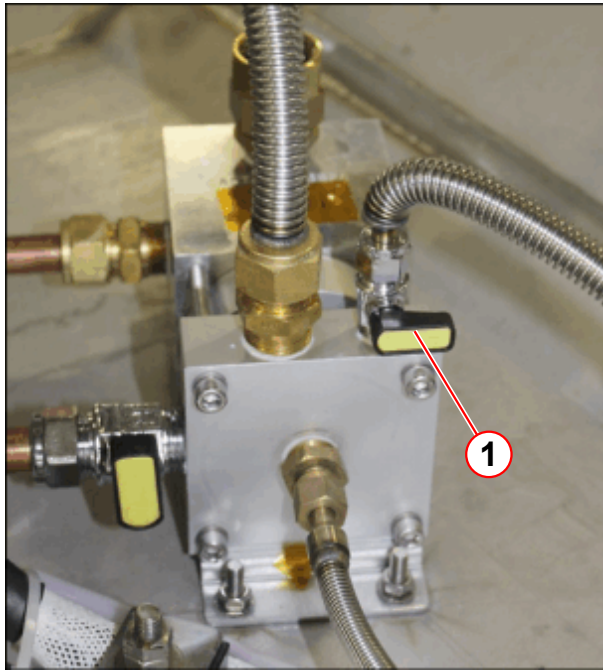

(1) Bypass valve


- The cold head vent line is to remain connected, with the hand-valve **CLOSED** during normal system operation.

Fig. 72: Leak checking the cold head seals

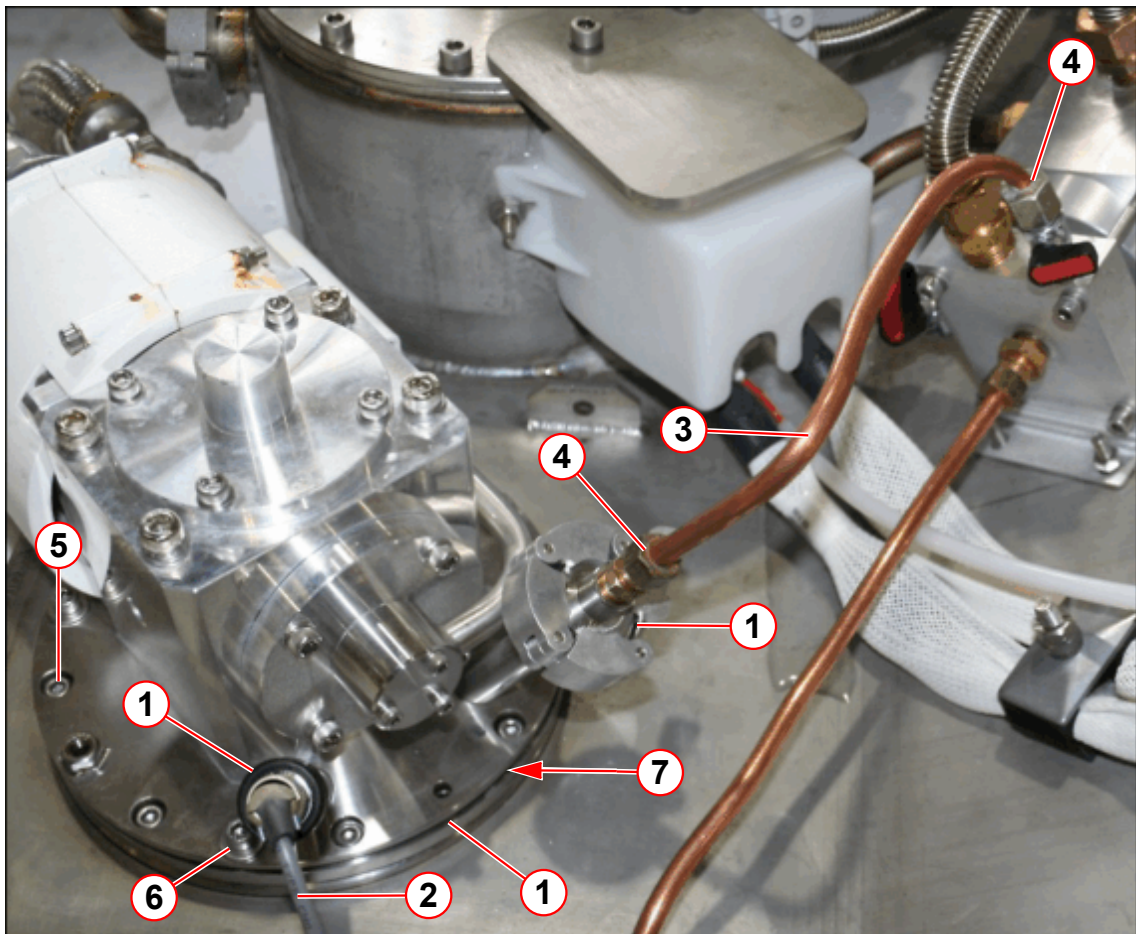


- (1) Cold head seals leak check points
- (2) Temperature sensor cable
- (3) Cold head vent line
- (4) Cold head bypass line connection
- (5) 8 off M6 x 16 mm O-ring clamp ring screws, torque to = 8Nm
- (6) 4 off M6 x 25 mm cold head flange screws, torque to = 4Nm
- (7) O-ring seal

Fig. 73: Cold head 10 pin seal

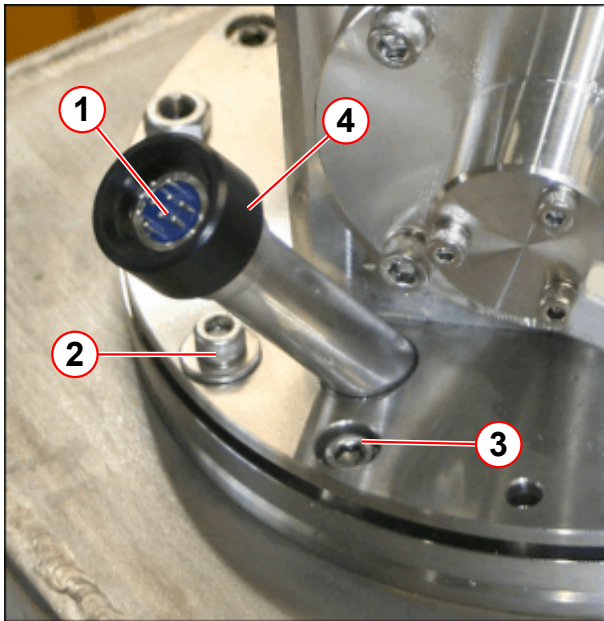

- (1) Leak check point
- (2) 4 off M6 x 25 mm cold head flange screws, torque to 4Nm
- (3) 8 off M6 x 16 mm skt hd screws O-ring clamp plate, torque to 8Nm
- (4) Black sealing cap

Fig. 74: Protective heat shrink &amp; warning label

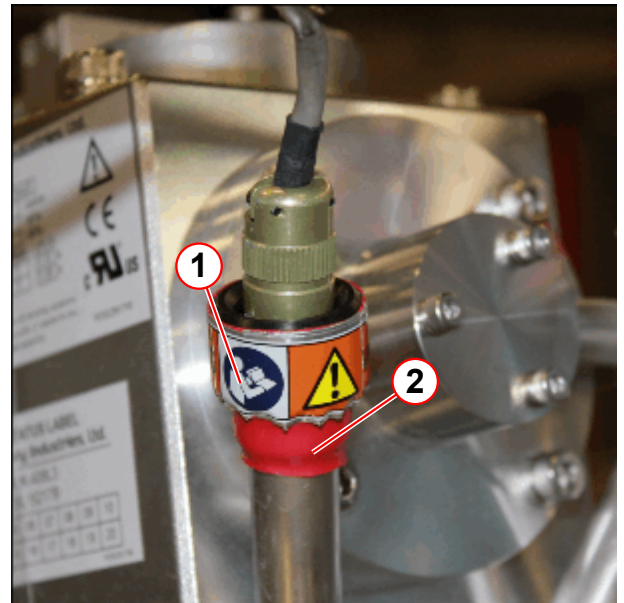

- (1) Warning label
- (2) Heat shrink



If a leak is found with the magnet at field, it is permissible to re-torque the 4 off cold head flange screws to a higher setting of 5 Nm.

If the leak persists, send a report to your local Uptime Service Center for repair by a specialist.

### 3.4.6.2 Leak Checking the Temperature Sensors



#### WARNING

Risk of cold burns. High helium gas pressure.

Risk of personal injury may occur.

- » Do not unscrew or remove the black temperature sensor sealing cap.
- » Remove the temperature sensor cable with a straight pull away from the sealing cap.




Cold heads fitted from 2015 have a protective heat shrink and a warning label fitted around the temperature sensor connector.

- **DO NOT** turn or try to remove the 'black' sealing cap.
- Remove the temperature sensor cable with a straight pull away from the black sealing cap.
- Leak check the 10 pin seal connector.




- Immediately refit the cable if no leak is found.

Fig. 75: Temperature sensor

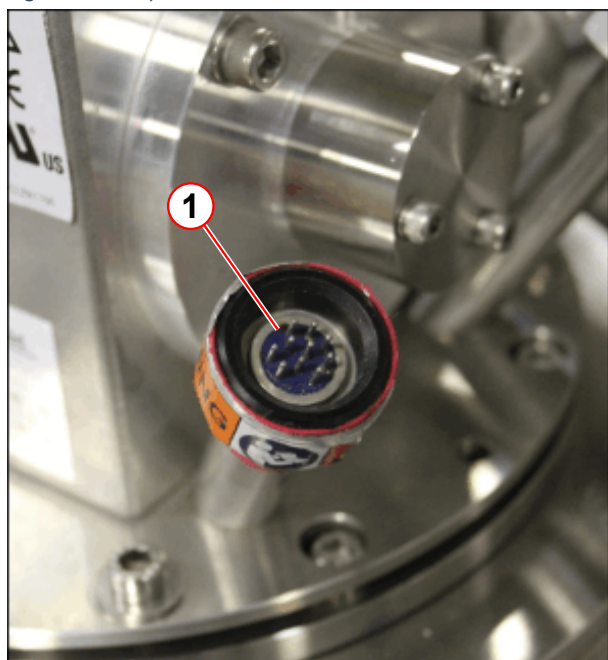

(1) 10-pin seal connector


Fig. 76: Protective heat shrink & warning label

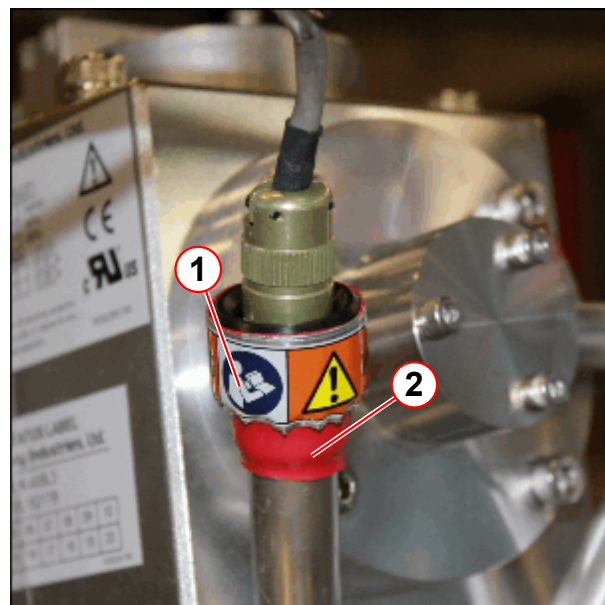

(1) Warning label  
(2) Heat shrink

### 3.4.6.3 Syphon Port

Fig. 77: Syphon blanking plug

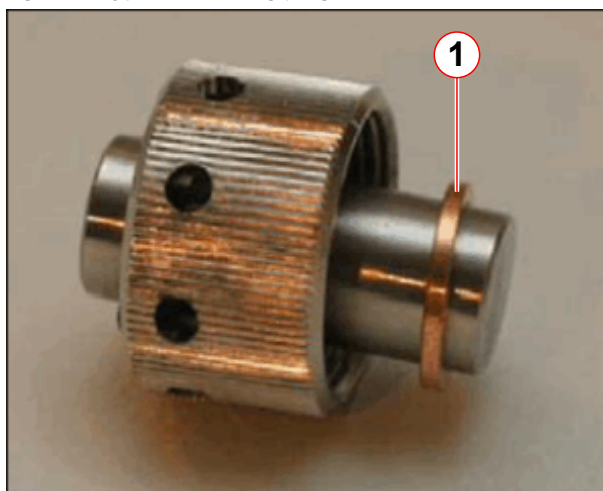


(1) Copper seal

- After all helium fills warm up the syphon port. and fit a new, copper seal, (material number **03108094**).
- Tighten the ring nut firmly. Do not over- tighten the syphon plug.



Fig. 78: Syphon fill port

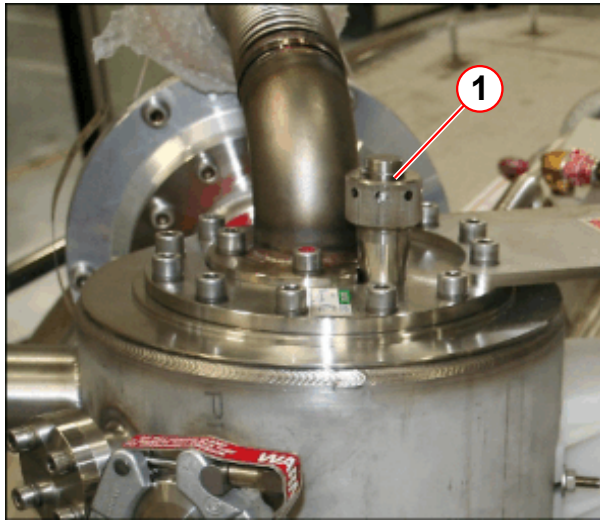

(1) Syphon port plug fitted


- Warm up and dry around the syphon port.
- Tighten the ring nut again. Do not over-tighten the syphon plug.
- Leak check. Repeat warm up of syphon plug and tightening of ring nut.
- If a leak persists, replace the copper seal.

#### 3.4.6.4 Pressure Gauge Assembly

Leak check the pressure switch transducer assembly.

- Leak check around the fixings of the vent line to the pressure gauge.

Fig. 79: Leak check the pressure line

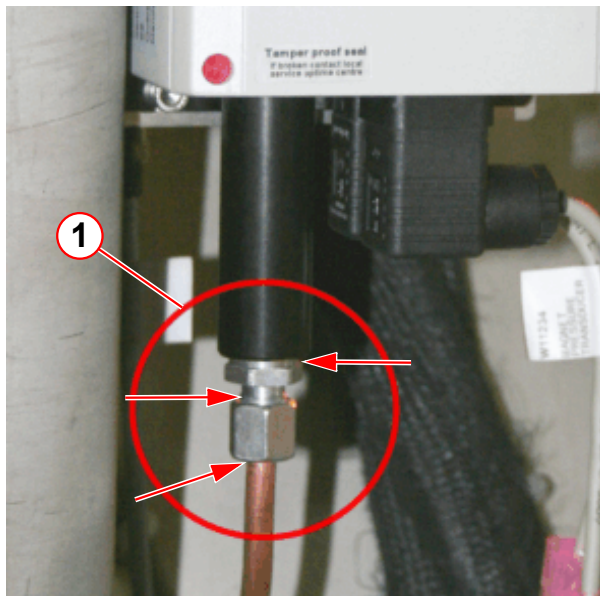

(1) Leak check points

- For leaks found at pipe couplings, first try tightening the compression nut and high load seal. **BEWARE OF OVER-TIGHTENING!**
- If a leak persists the compression seal must be replaced.

#### 3.4.6.5 FCL Lead Connector

The FCL lead connector passes through the turret in a ceramic interface. It is made leak tight by a high load seal. The FCL ceramic connector and seal should only be checked if no other leaks have been found.

- A broken ceramic seal will be seen as a gross leak.

- Report a damaged ceramic seal to the Magnet HSC. The FCL turret connector and internal wiring may need replacing.


## NOTICE

### Air ingress

The ceramic seal is a fragile component

- » Take care when removing the FCL cover or working next to the ceramic seal.

Fig. 80: FCL cover and location

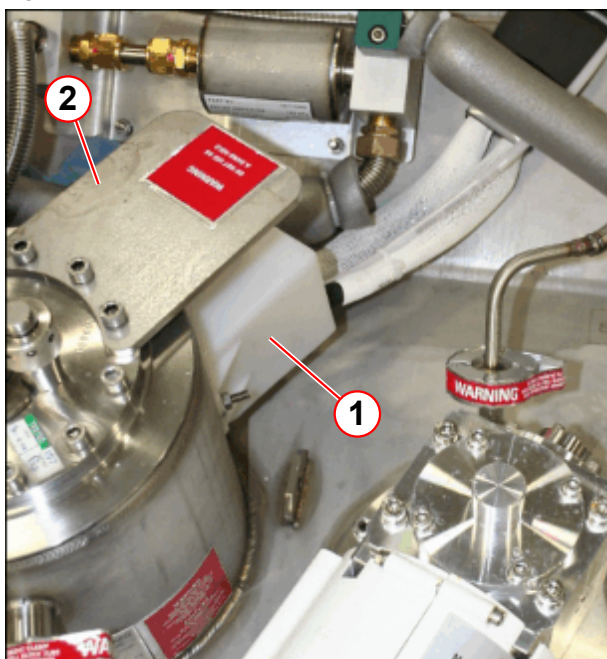


- (1) FCL cover
- (2) FCL protection plate

The FCL connector is located under the plastic protection cover.

- Only remove the cover if a leak is detected around this area.

### 3.4.6.6 16 / 17 psi (A) Vent Valve Assembly

Fig. 81: Vent Valve - High pressure points

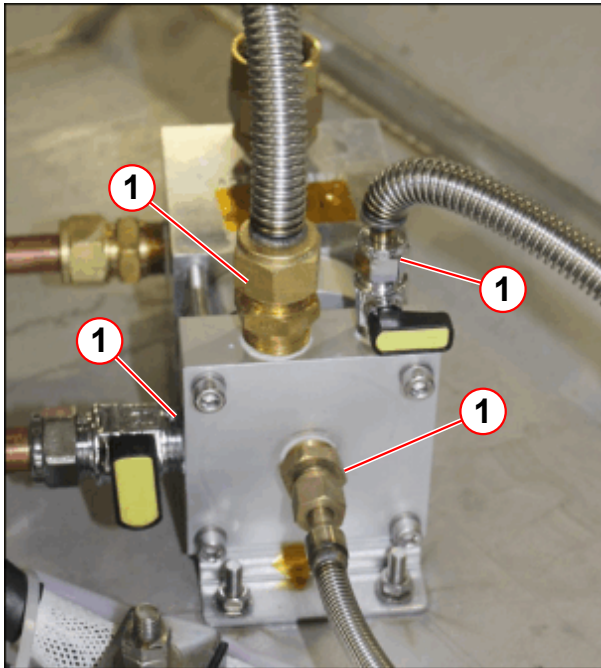


(1) Leak check locations

- Leak check all high-pressure joints on the magnet vent valve assembly (Pos.1).
- For leaks found at pipe couplings, first try tightening the clamp or securing screws. **BEWARE OF OVERTIGHTENING!**

### 3.4.6.7 Turret and Venting

The figure below shows the high pressure leak check turret locations.

Fig. 82: Turret leak check - high pressure

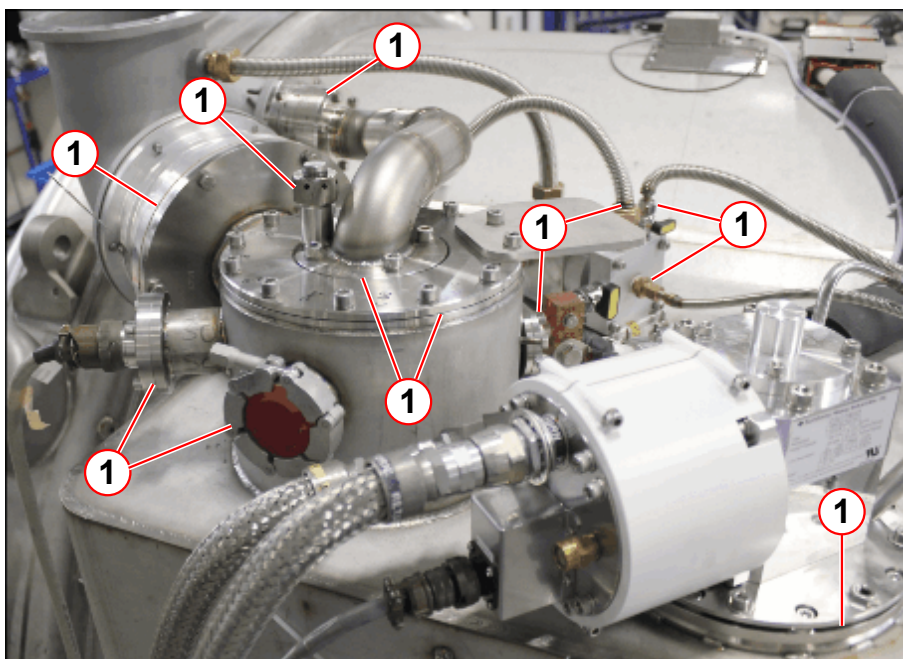

(1) Leak check points

### 3.4.7 Internal Magnet Leak Checks

Leaks can occur through the magnet vent valve or burst disks:

- Auxiliary bursting disk.
- Quench valve burst disk assembly.
- 17 psi (A) valve assembly.



#### CAUTION

Internal leak checks.

Risk of air ingress causing internal icing of turret.

- » Replace burst disk or valve immediately if found to be broken or leaking.



The magnet must be at **ZERO FIELD** before leak checking a burst disk.

**WARNING****Risk of air ingress**

If a burst disk is broken there is a high risk of air ingress and ice.

- » **IMMEDIATELY** fit the replacement burst disk from the spare parts kit.
- » Order a replacement burst disk.

**3.4.7.1 Auxiliary Bursting Disk**


The auxiliary burst disk should not break during a magnet quench.

- If necessary warm up the burst disk and venting with a heat gun.
  - Activate the pressure heater.
  - Allow the system to pressurize.

Leak check externally around the high pressure seal (pos.1).

- If there is an internal leak or breakage the magnet pressure will be zero.
- The magnet must be at zero field before further checks and replacement procedures can be performed.

Fig. 83: Auxiliary burst disk assembly

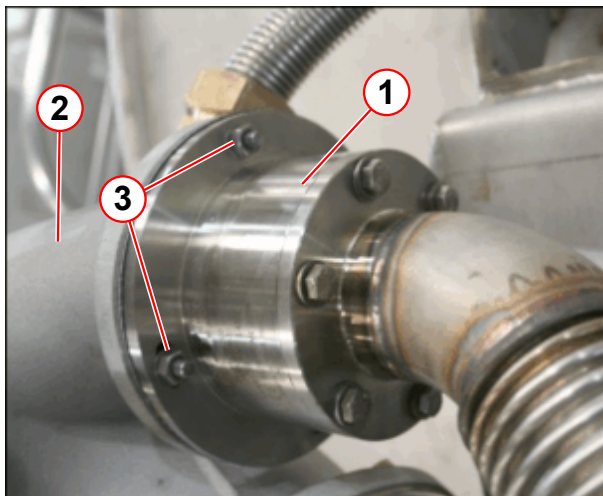

- (1) Auxiliary burst disk
- (2) Elbow
- (3) M6 screws, washers and nuts (4 off)

To check for any leak path through the auxiliary metal bursting disk:


- Disconnect the auxiliary burst disk from the elbow, on the low pressure side.
  - Undo the 4 x M6 screws, washers and nuts (pos.3).
- Carefully move the auxiliary vent pipe away from the elbow.



Do not touch the metal membrane; it is very fragile.



Fig. 84: Auxiliary burst disk

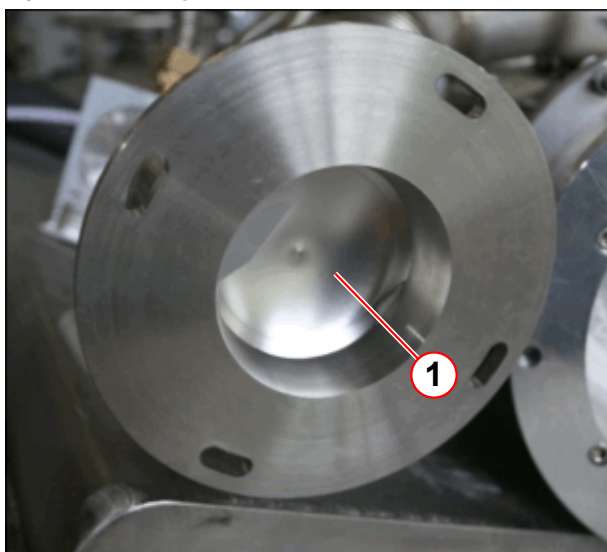

(1) Burst disk

- Visually inspect the internal membrane of the burst disk for damage.
- Visually inspect the O-ring on the low pressure side.
- Replace the auxiliary metal burst disk and high load C-seal if a leak is detected.

#### 3.4.7.2 Quench Valve Burst Disk Assembly




The magnet must be RUN TO ZERO before checking the quench valve.

If the magnet has a low gauge pressure and will not pressurize, it is possible that the quench valve has an internal leak.

The quench valve burst disk and spring plate must be leak checked at installation & Start Up and Post quench.

- Remove the quench valve elbow assembly from the turret to access the burst disk.


Fig. 85: Vent valve outlet line

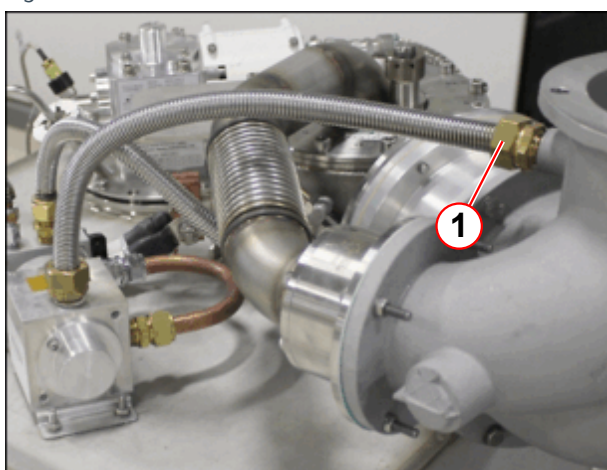

(1) Vent line elbow connection

- Disconnect the vent valve outlet line from the quench line elbow.


Fig. 86: Vent line connections

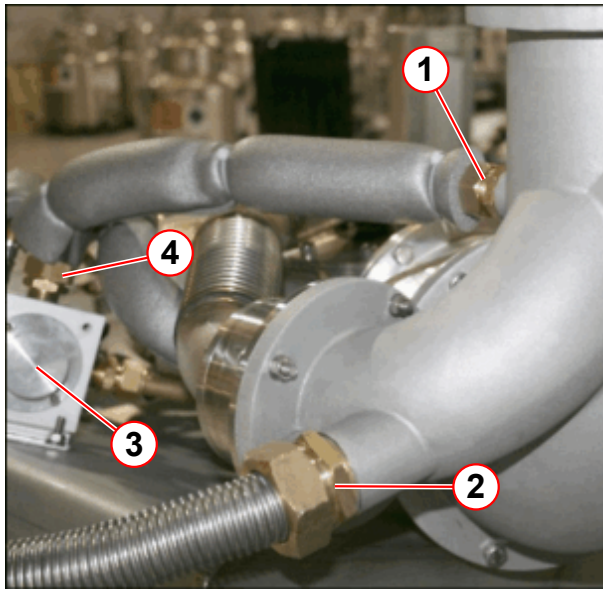

- (1) Elbow outlet line
- (2) OVC outlet line
- (3) Vent valve assembly
- (4) Vent valve outlet line

- Disconnect the OVC outlet line from the elbow (**pos.2, if fitted**).

Fig. 87: Vent valve connection

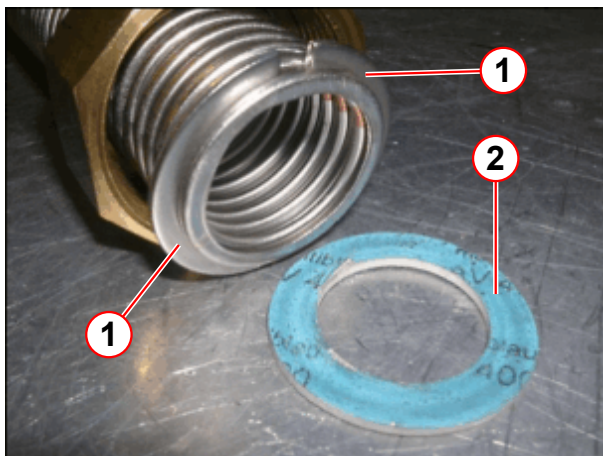


- (1) Split rings
- (2) Gasket

Take care not to lose the split rings during this process!

- Ensure the split rings do not fall out of position.
- Be careful to retain the gasket when disconnecting the vent line.
  - The gasket is not retained by the split rings and nut.



Fig. 88: Quench valve &amp; burst disk

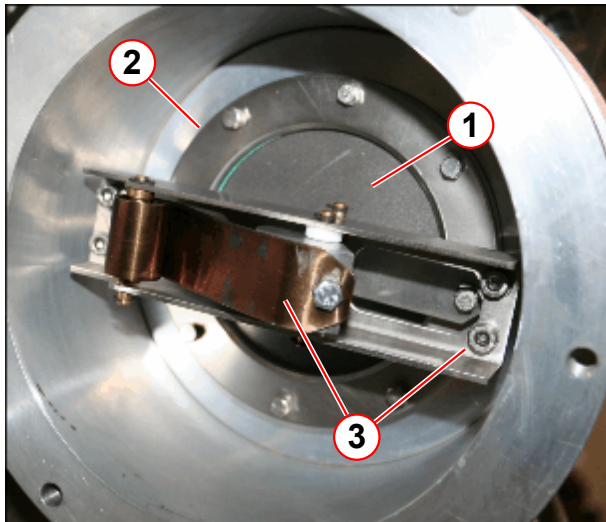


- (1) Graphite burst disk
- (2) Quench valve O-ring (high pressure)
- (3) Quench valve spring & screws

Perform a leak check of:

- The graphite burst disk.
- Around the O-ring seal to the turret (pos.2).
- Check the condition of the spring and retaining screws.

Fig. 89: Quench valve burst disk and O-ring

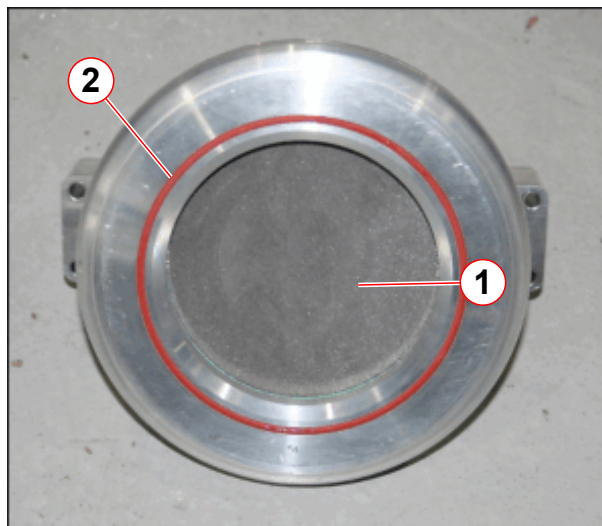

- (1) Burst disk
- (2) O-ring

If a leak is seen around the burst disk or spring plate the magnet must be depressurized and the burst disk assembly removed.


- » Details of removing the quench valve is in the Replacement procedures.
- Inspect the burst disk.
  - » If the burst disk has ruptured or is cracked, it must be replaced.
- Inspect the burst disk O-ring.
  - » If the O-ring is split or damaged it must be replaced.
  - » O-ring material number **10122816**.

### 3.4.7.3 16 / 17 psi (A) Vent Valve Assembly

If the magnet has a low gauge pressure and will not pressurize, it is possible that the vent valve assembly has a gross internal leak.

- Check the magnet vent valve assembly for an internal leak.
- If a leak is suspected the magnet must be ramped to zero field.


Fig. 90: Vent valve outlet line

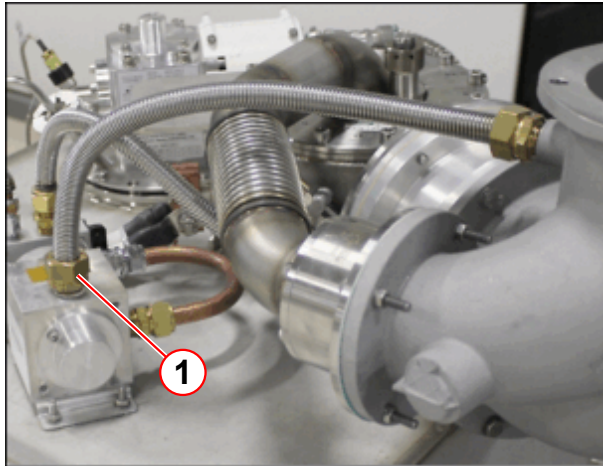

(1) Vent line connection

- Remove the low pressure vent line from the vent valve.

Fig. 91: Vent valve connection

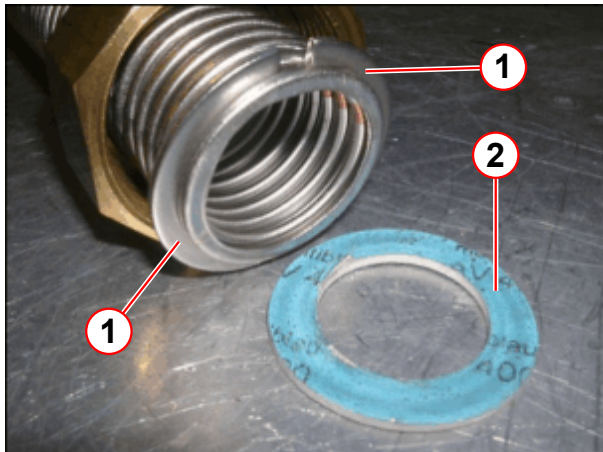


(1) Split rings  
(2) Gasket

Take care not to lose the split rings during this process!


- Ensure the split rings do not fall out of position.
- Be careful to retain the gasket when disconnecting the vent line from the elbow.
  - The gasket is not retained by the split rings and nut.

Fig. 92: Glove fitment

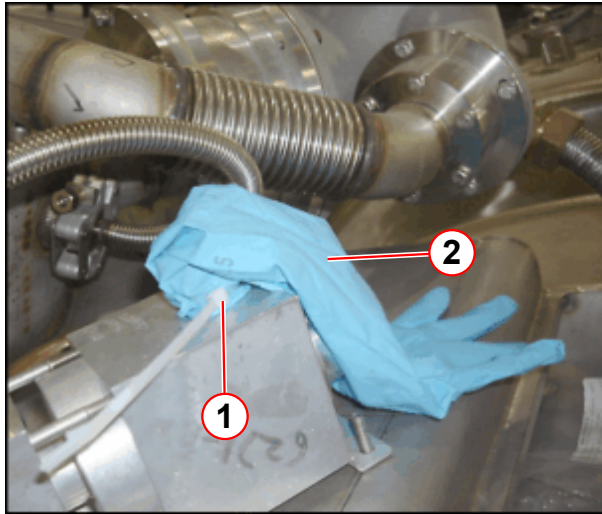

- (1) Surgical glove
- (2) Cable tie

Fig. 93: Gross leak through the vent valve assembly

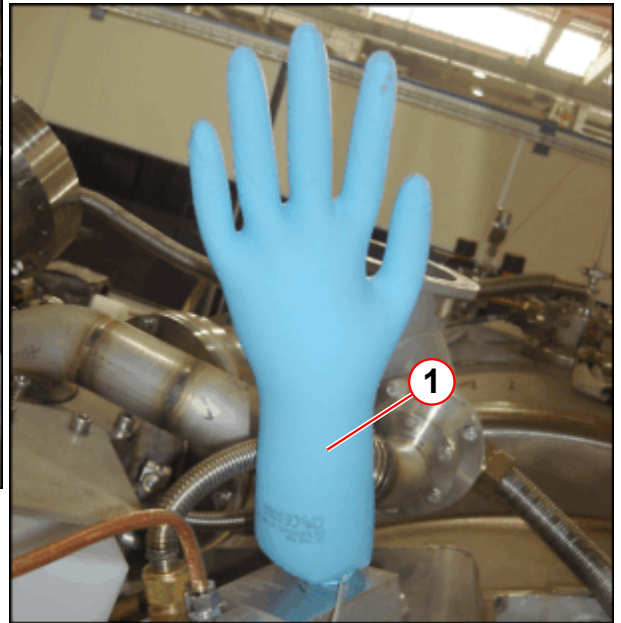

- (1) Glove inflated

- Fit a surgical glove (or small airtight plastic bag) over the end of the open low pressure vent.
- Seal the glove with a cable tie.
  - » If the glove expands rapidly with the magnet pressure lower than 15.6 psi (A) (0.3 psi (G)), this indicates a gross leak.
  - » If the valve is not sealing, replace the bypass valve assembly.


### 3.4.8 Low Pressure Turret and Venting Components



Low pressure venting leaks can only be checked visually during depressurization, filling or a magnet quench.

- Check for cold gas escaping around low pressure joints during:
  - depressurization, helium filling or a magnet quench.
- Check that all sealing nuts and flange screws on the venting are tight.
  - Warm venting and replace low pressure gasket seals and O-rings if found damaged.


Fig. 94: Low pressure leak checks

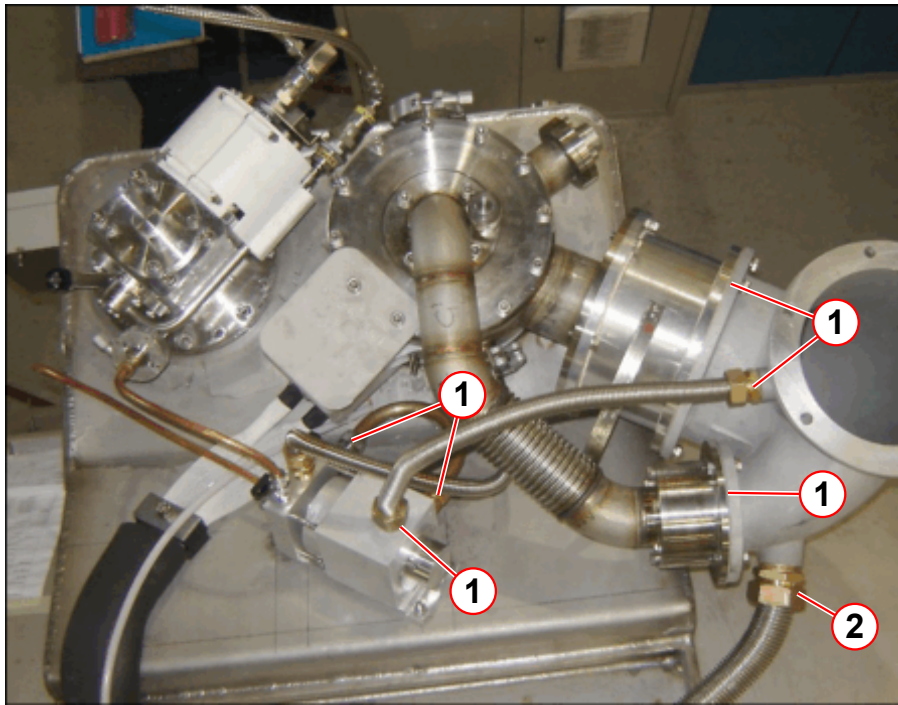

- (1) Low pressure joints
- (2) OVC burst disk port (if fitted)

Low pressure venting leak check locations with no OVC (DOP port) vent line connection to the elbow.


Fig. 95: Turret and venting leak check

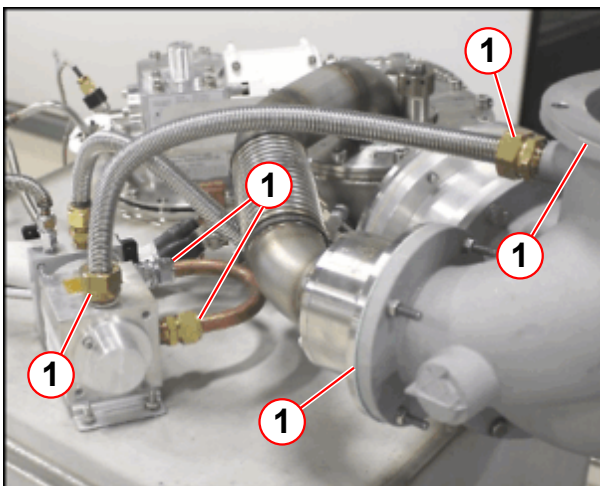

- (1) Low pressure points



Fig. 96: Horizontal service port

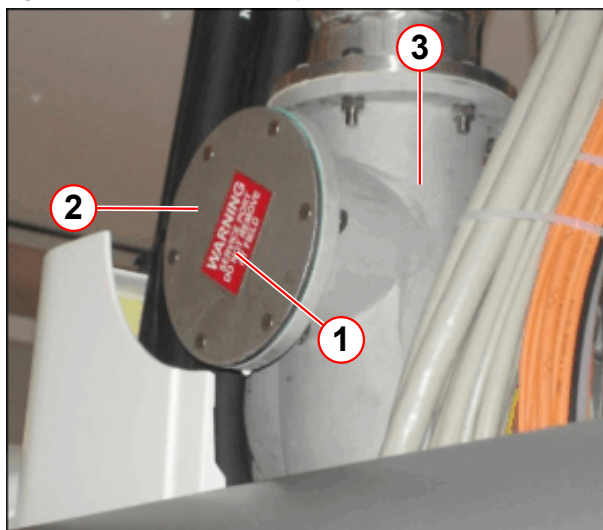

- (1) Warning label
- (2) Port (blanked)
- (3) Quench elbow

#### OR98 only - Horizontal Service Port

- This port cover must NOT be removed.
- **Do not** use as an inspection cover to check the quench valve burst disk.
- Check that the low pressure gasket seal is leak tight.

### 3.4.9 Final Work Steps

- After the leak check, ensure that the leak detector is working correctly:
  - » Check the equipment function with a known helium gas source.
- Ensure that all parts of the venting system have been connected and fitted.
- Check that the magnet has returned to the normal operating pressure.
- Refit all system covers.
- Order replacement seals if spares have been used.

## 4.1 OR98/99/103 Ramping



Only Siemens CSEs trained in ramping magnets are allowed to perform the magnet ramping procedure.



Before ramping check that the quench line is correctly installed from the magnet to the outlet.



### WARNING

#### Magnet safety checks

Failure to observe the following safety warning may result in physical injury and equipment damage.

- » Ramping is not allowed if the visually inspection of the quench line was not performed successfully.
- » Before initially ramping up the magnet, the 'Testing the ERDU circuit' has to be performed.
- » Procedure described in the safety related test can be used if ERDU test load is available.



### WARNING

#### Strong magnetic field

Failure to observe the following safety warning may result in physical injury and equipment damage.

- » Before ramping up the magnet, ensure that no ferrous tooling or equipment are present in the imaging suite.
- » Ensure that the warning signs identifying the controlled access area (0.5 mT zone) have been posted at the entrance to the RF room as well as all other areas leading to the 0.5 mT zone.
- » The sign must be posted in a way that it is visible even when the door to the controlled area is opened.



### CAUTION

Risk of asphyxiation and hypothermia, if helium gas is vented into the magnet room

Failure to observe the following safety warning may result in dizziness or loss of consciousness in personnel.

- » Quench valve must be connected to the Quench pipe.



Initial installation ramp helium level must be:

- ≥ 60 % for the OR98
- ≥ 65 % for the OR99 and the OR103.

Ensure a minimum of:

- OR99/OR103 = 2000 litres of helium on standby.



During normal Service ramping the minimum helium level for ramping up must be:

- > 40 % for the OR98 and the OR103
- > 60 % for the OR99.



The full functionality of the MRC (MRAWP) Host computer is not available during ramping procedure.

Do not re-boot the system during ramping.



Wait at least one hour after completing a previous ramp before starting the next ramp, (including lead checks and aborted ramps).



Wait at least two hours after helium filling before ramping the magnet.



Check that the venting bypass valve is in the CLOSED position.

#### Always

- Always restrict access to the magnet room when ramping the magnet. Use flashing lights and signs.
- Always inspect all equipment and tools for defects prior to use.
- Always inspect the Anderson ramp lead connections at the filter plate.
- Clean the area of the magnet room, dispose off empty cartons or other packing materials
- Always inspect magnet room for the presence of non identified objects. Before ramping, the nature of those objects must be clarified.



Prepare a minimum of 4 pictures of the room **prior to first ramping**.

These pictures must be colored and should have 2 to 4 megapixel resolution. They must show the complete room, all 4 corners and walls. In case there are un-removable cupboards installed, their doors must be open to prove that they are empty.

Later those pictures are to be attached to the installation report.



**Never**

- Never use defective or damaged equipment.
- Never allow unauthorized people to stand too close to the magnet during ramping.
- Never block emergency exits.
- Never smoke or use open flames in the magnet room.
- Never ramp a magnet when you are unsure about magnetic parts inside the magnet room. Those parts could be hidden in e.g. supposed empty boxes or in closed cupboards or sliding of a desk.

#### 4.1.1 Tools required

**WARNING****Strong magnetic field**

Failure to observe the following safety warning may result in physical injury and equipment damage.


- » Only non-magnetic equipment may be taken into the magnet room.

- Magnet Power Supply Unit (MPSU 3600); (order number: **8395613**). This includes:-
  - 3600 MPSU Power lead.
  - CAN Open lead.
  - RS232 lead.
- Heat gun.
- MRC (MRAWP) Host or Service PC.

### 4.1.2 Preparation

- Move the 3600 MPSU into the equipment room close to the RF filter plate.

Fig. 97: MPSU

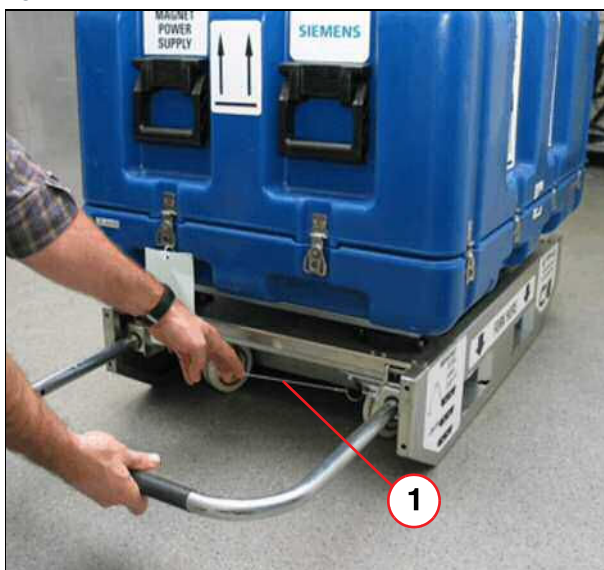

(1) Wire locking strap

- Lower the trolley handle into the horizontal position. This will raise the wheels and allow the trolley to rest on the two skids.
- Pull the wire strap to unlock the handle, and slide the handle into the recess under the trolley.
- Remove the shipping case lid by unclipping the 12 spring clips.

Fig. 98: MPSU



(1) MPSU LID

- Open the hinged lid of the 3600 MPSU by pulling the two silver knobs and turning them through a 1/4 turn.
- Raise the lid.

### 4.1.3 Installation

Fig. 99: 3600 MPSU Power connections



- (1) X1 (MPSU power)
- (2) CB1 (Mains circuit breaker)

- Ensure the mains circuit breaker **CB1** is in the **OFF** position on the 3600 MPSU.

Fig. 100: Adapter plate (view from Equipment Room)



- Unscrew the securing screws and remove the **adapter plate** from the **RF filter plate**.
- This will allow access to the magnet ramp leads in the **magnet room**.



Do NOT connect the 3600 MPSU ramp leads to the MPSU until AFTER the Transport test has been completed. See (→ Transport Test Procedure / Page 117).



The mandatory transport test was introduced from 3600 MPSU software revision level 17.

This test is only required the FIRST TIME time the 3600 MPSU is used on site.

Tab. 8 MPSU lead connections

| Lead              | From               | To                                                                                                                                           | Length |
|-------------------|--------------------|----------------------------------------------------------------------------------------------------------------------------------------------|--------|
| Magnet Ramp leads | MPSU X3<br>MPSU X4 | Filter plate (grey)<br>Filter plate (blue)                                                                                                   | 5m     |
| CAN OPEN lead     | MPSU X7            | EPC X125                                                                                                                                     | 5m     |
| RS232 lead        | MPSU X6            | MRC (MRAWP) Host, or Laptop<br><b>W10220 / X12 / MPSU</b><br>Cable W10220 to <b>X128</b> , disconnect cable at <b>X12</b> on top of the EPC. | 5m     |
| MPSU Power lead   | MPSU X1            | EPC <b>X130</b>                                                                                                                              | 10m    |

Fig. 101: OR98, OR99 and OR103 Ramping configuration


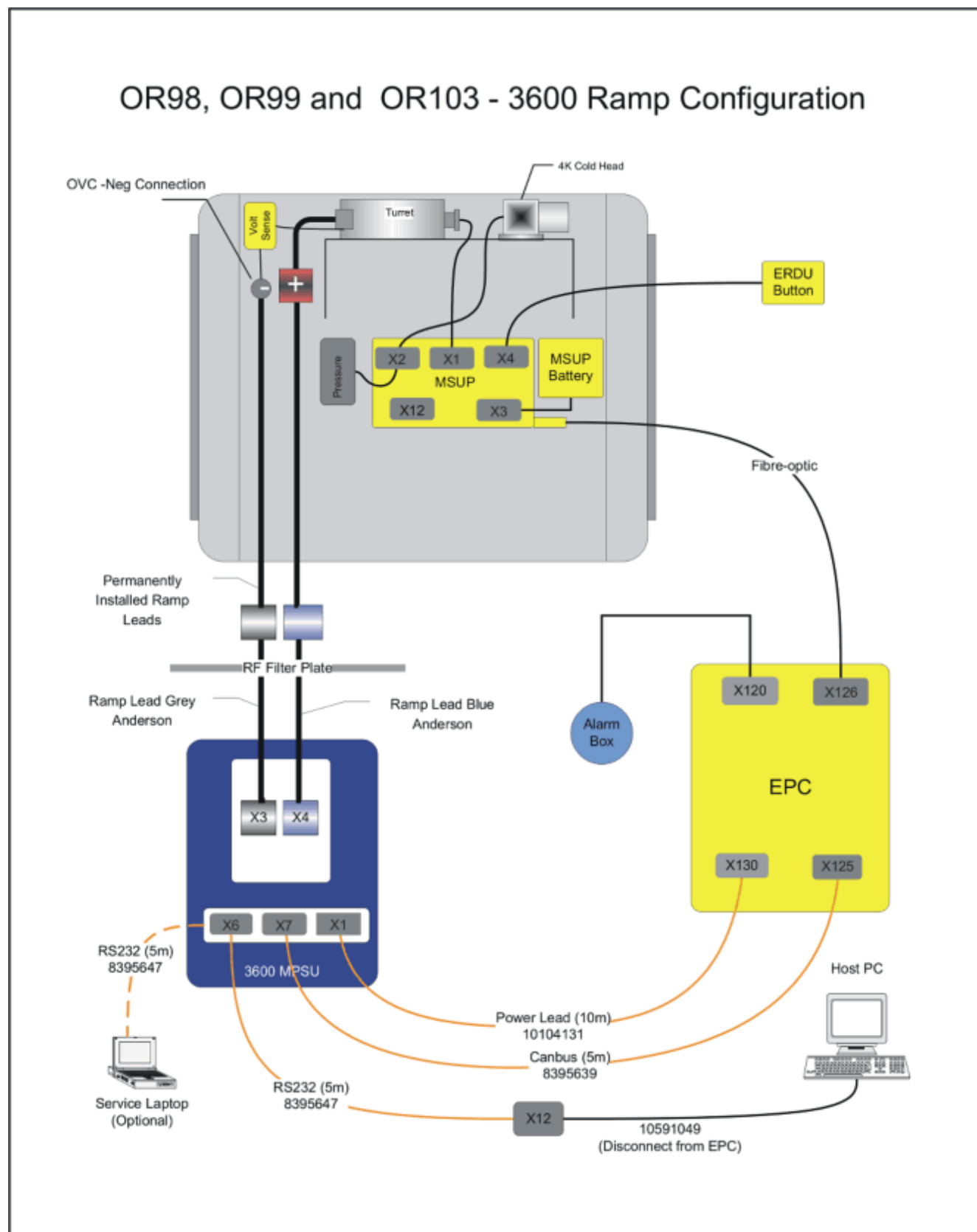


Fig. 102: 3600 MPSU Connector and indicator panel




- (1) ON indicator
- (2) FAULT indicator
- (3) SELF TEST indicator
- (4) CONTACTOR OPEN indicator
- (5) CURRENT IN LEADS indicator
- (6) SWITCH HEATER ON indicator
- (7) X7 CAN OPEN connector
- (8) X6 RS232 connector

Fig. 103: 3600 MPSU leads

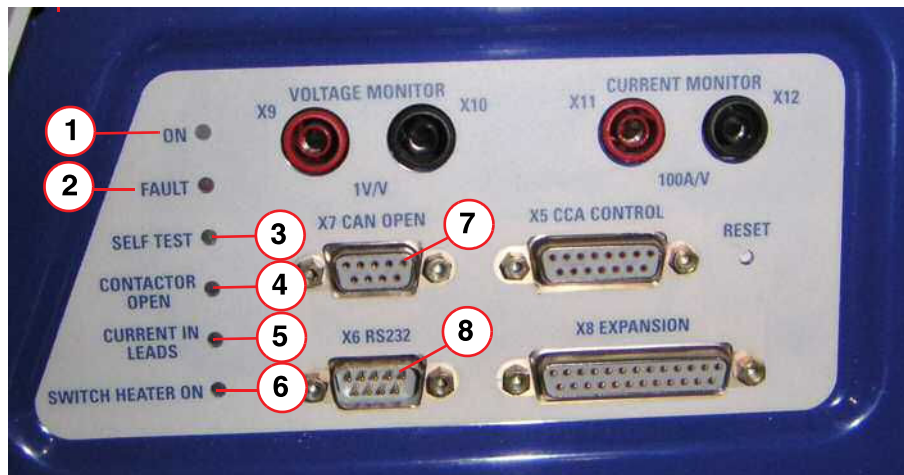

- (1) CAN OPEN lead
- (2) RS232 lead
- (3) MPSU power lead

Fig. 104: EPC connections

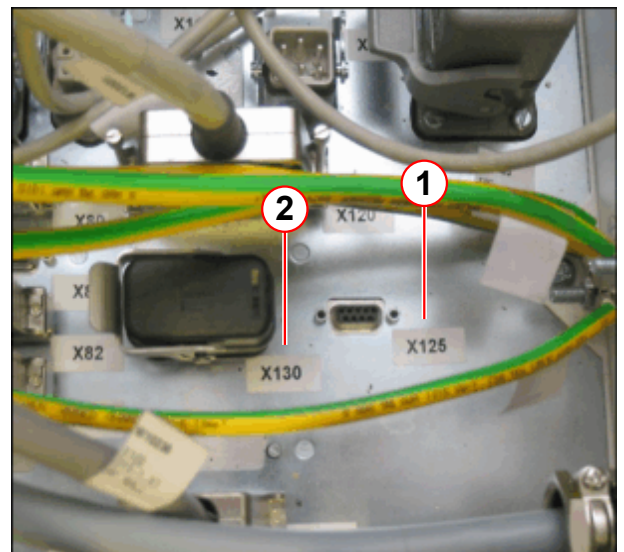


- (1) X125
- (2) X130

- Connect the **CAN OPEN**, **RS232** and power leads. (→ Fig. 103 Page 99)
- Ensure the system is switched **ON** on the **alarm box**.



Fig. 105: MPSU Circuit breaker F53

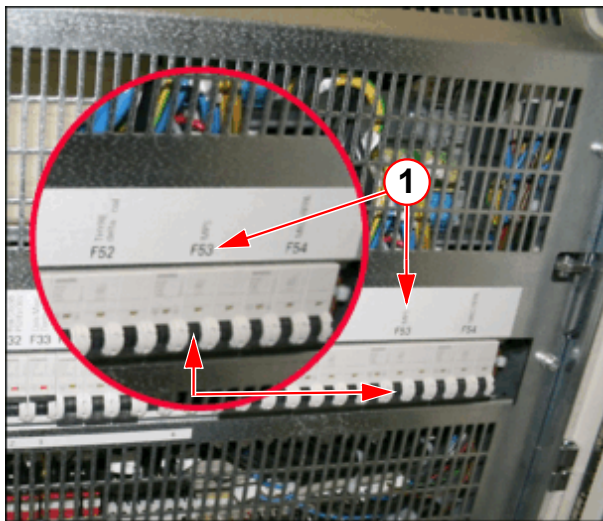

(1) Circuit breaker F53

- Turn ON **F53** in the **EPC**.

#### 4.1.4 Set up of the MPSU




To start ramping, the **MRC Host** or **Service laptop** must be connected to **X6** of the 3600 MPSU (RS232 interface) (→ 8/ Fig. 102 Page 99).

##### 4.1.4.1 Service Laptop connection

Power to a Service laptop is supplied at **X2** on the 3600 Power Supply.

Fig. 106: 3600 MPSU Fuses and laptop power supply

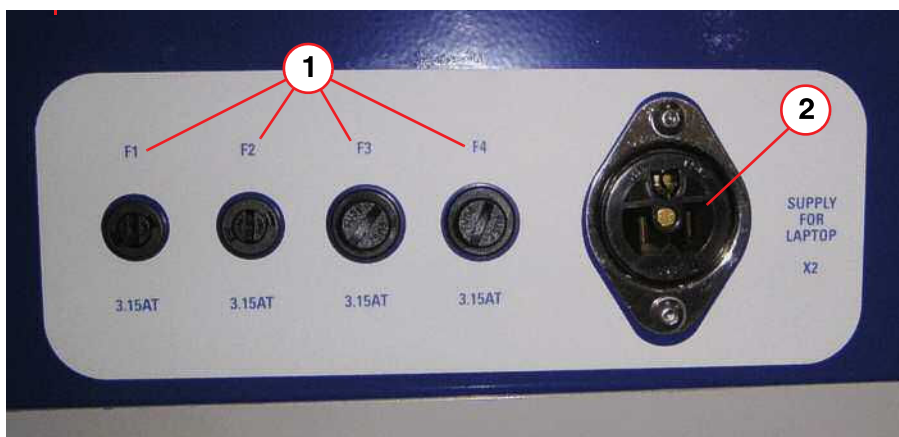

- (1) Fuses F1 - F4
- (2) X2 Supply for Laptop computer

Establish communication using a **Terminal emulation software** (e.g. PuTTY) (→ Using Putty / Page 102)



## 4.1.4.2 MRC Host connection



The following work step is only necessary with SW N4\_VD11A.

With SW N4\_VD11D, the ttyrecorder stops automatically before the terminal emulator starts, and re-starts it when the terminal emulator is closed.

- Hold down CTRL/ESC together
  - » select **programs - advanced user - OK**
- Type in password
  - » meduser1
- Select OK
- Select Start - Run
  - » The menu appears
  - » Type in CMD
- Select OK


Before running **Terminal emulation software**, the Host process **stop\_ttyrecorder.bat** must be run. This will release the serial port to enable the host to communicate with the MPSU. This process is automatically re-started after a reboot.

- In the DOS window type **stop\_ttyrecorder.bat** <enter>

## 4.1.4.3 Configuring the Host

- Open the **Magnet & Cooling** service platform.
- Select **Terminal** from the **toolbar**.

Fig. 107: Open the Terminal via Magnet&Cooling

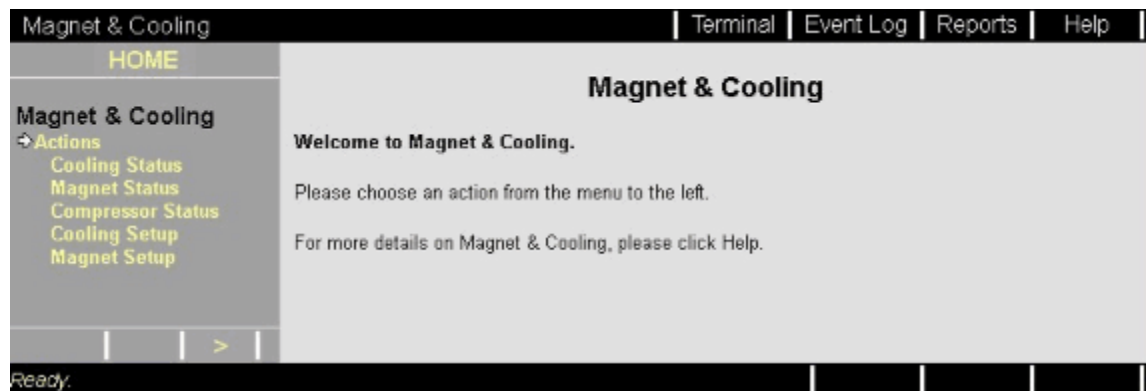

- **Terminal** option: the program will **open**. No configuration is required.
- **Terminal** emulator option: the program will open, configuration is required


## 4.1.4.4 Using Putty

## Setting up PuTTY

You can use software called PuTTY to connect to MSUP's (Magnet Supervisorys) and MPSU's (Magnet Power Supplies). This software has to be setup with the following settings in order to work.


Ensure you have PuTTY version 0.60 or newer software installed on your PC.

Fig. 108: PuTTY Configuration - Session

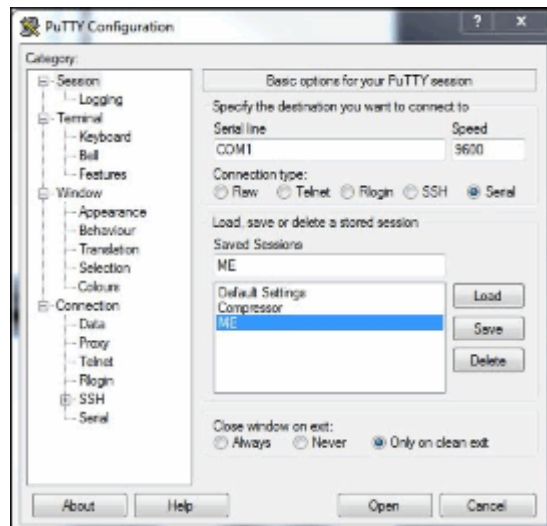


- Under Session tab.
- Set the 'Connection type:' to Serial.
- Set the 'Serial line' to the correct COM port being used. In this case COM1. Only state this if the com port that is being used. Change as required.
- Set the 'Speed' to 9600.
- Set the 'Close window on exit' to Only on clean exit.

Fig. 109: PuTTY Configuration - Logging

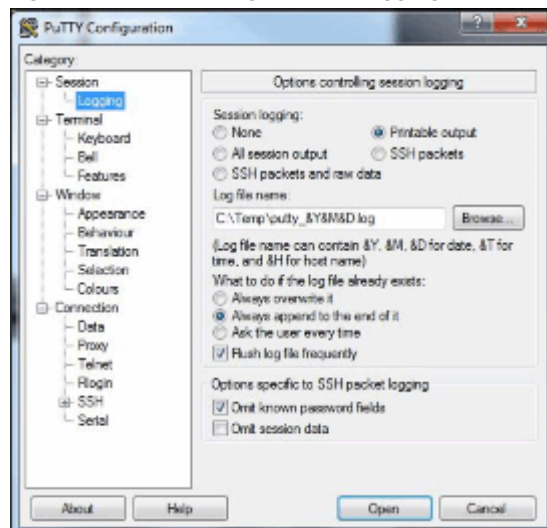


- Under Logging tab.
- Set the 'Session logging' to Printable output.
- Set the 'Log file name:' to C:\Temp\putty\_&Y&M&D.log.
- Set the 'What to do if the log file already exists:' to Always append to the end of it and tick Flush log file frequently.
- Tick the 'Options specific to SSH packet logging' to Omit known password fields.

Fig. 110: PuTTY Configuration - Terminal

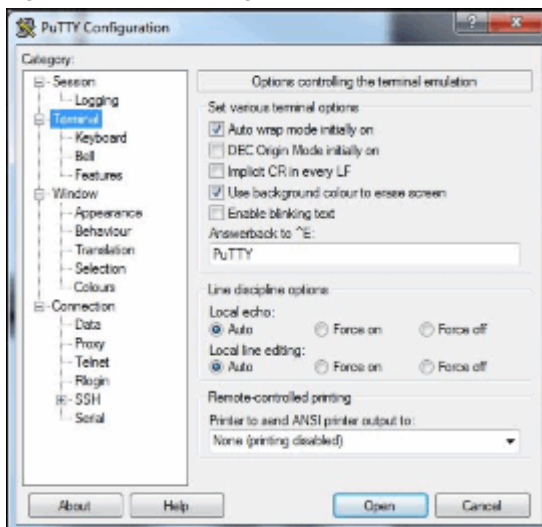


- Under Terminal tab.
- Under 'Set various terminal options' tick the following 'Auto wrap mode initially on' and 'Use background colour to erase screen'.
- Under the 'Answerback to ^E:' ensure it states 'PuTTY'.
- Under 'Local echo:' and 'Local line editing:' ensure 'Auto' is selected
- Under 'Printer to send ANSI printer output to:' select 'None (printing disabled)'.

Fig. 111: PuTTY Configuration - Keyboard

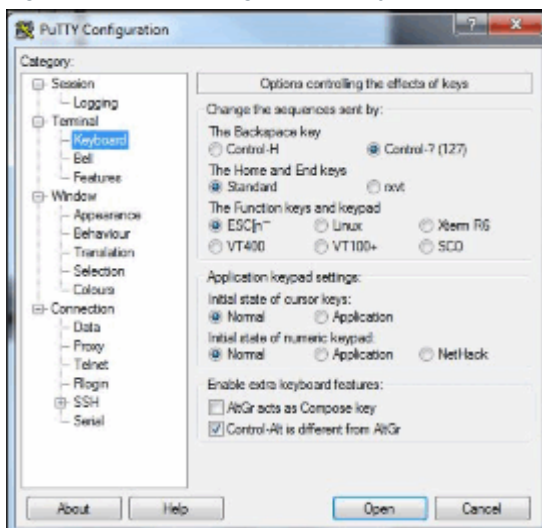


- Under the 'Keyboard' tab.
- Set 'The Backspace key' to 'Control-? (127)'.
- Set 'The Home and End keys' to 'Standard'.
- Set 'The Function keys and keyboard' to 'ESC{n~'.
- Set 'Initial state of cursor keys:' and 'Initial state of numeric keypad:' to 'Normal'.
- Tick 'Enable extra keyboard features' to 'Control-Alt is different is different from AltGr'.

Fig. 112: PuTTY Configuration - Bell

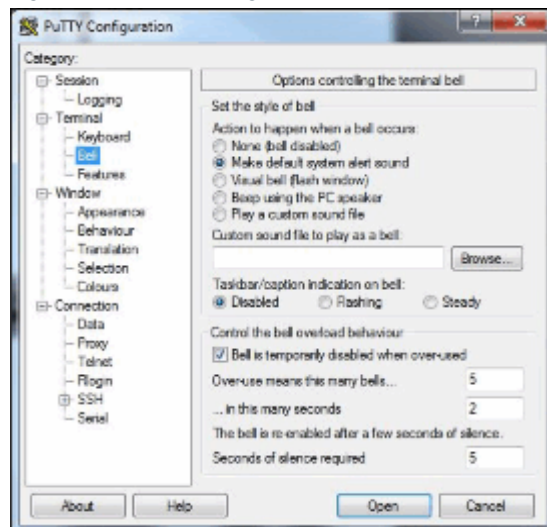


- Under the 'Bell' tab.
- Set 'Action to happen when a bell occurs' to 'Make default system alert sound'.
- Set 'Taskbar/caption indication on bell:' to 'Disabled',
- Tick 'Bell is temporarily disabled when over-used'.
- Set 'Over-use means this many bells...' to '5'.
- Set '...in this many seconds' to '2'.
- Set 'Seconds of silence required' to '5'.

Fig. 113: PuTTY Configuration - Features

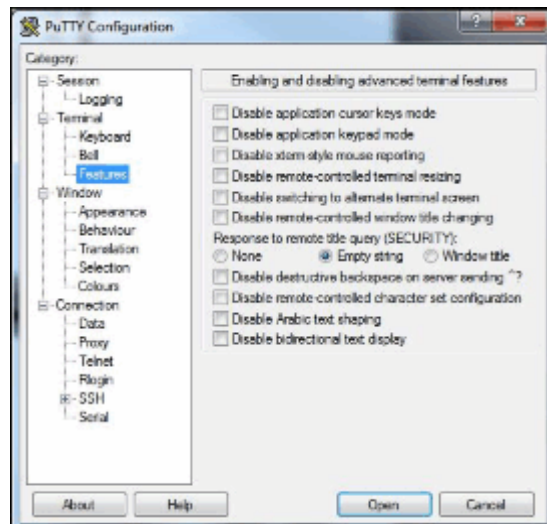


- Under the 'Features' Tab.
- Set 'Response to remote title query (SECURITY):' to 'Empty string'.
- Ensure nothing else is selected.

Fig. 114: PuTTY Configuration - Window

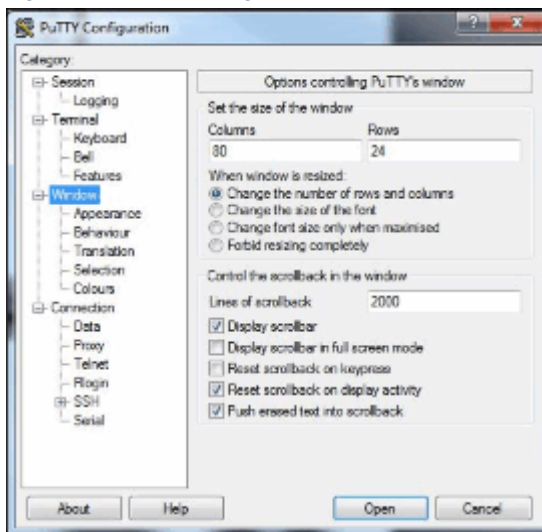


- Under the 'Window' Tab.
- Set the 'Columns' to '80' and the 'Rows' to '24'.
- Set 'When window is resized:' to 'Change the number of rows and columns'.
- Set the 'Lines of scrollbar' to '2000'.
- Ensure 'Display scrollbar', 'Reset scrollbar on display activity' and 'Push erased text into scrollbar' are all ticked.

Fig. 115: PuTTY Configuration - Appearance

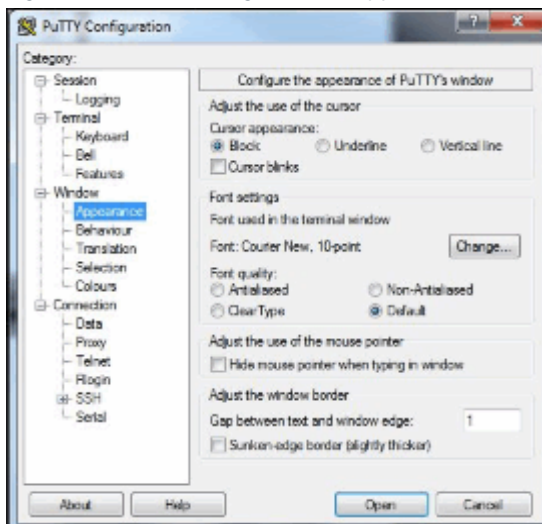


- Under 'Appearance' tab.
- Set 'Cursor appearance' to 'Block'.
- Set 'Font quality:' to 'Default'.
- Set 'Gap between text and window edge:' to '1'.

Fig. 116: PuTTY Configuration - Behaviour

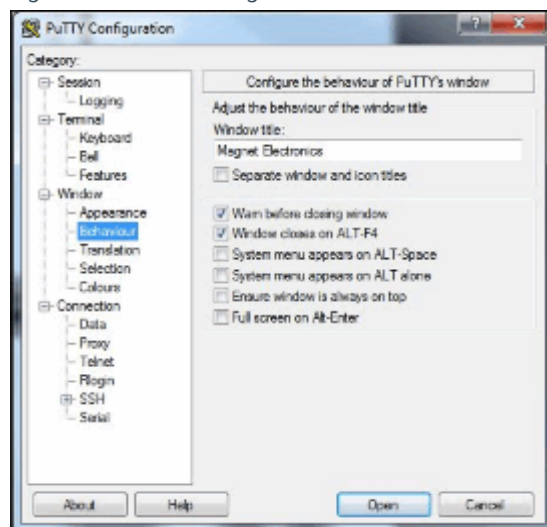

- Under 'Behavior' tab.
- Under 'Window title:' type 'Magnet Electronics'
- Ensure 'Warn before closing window', and 'Window closes on ALT-F4' are ticked.


Fig. 117: PuTTY Configuration - Translation

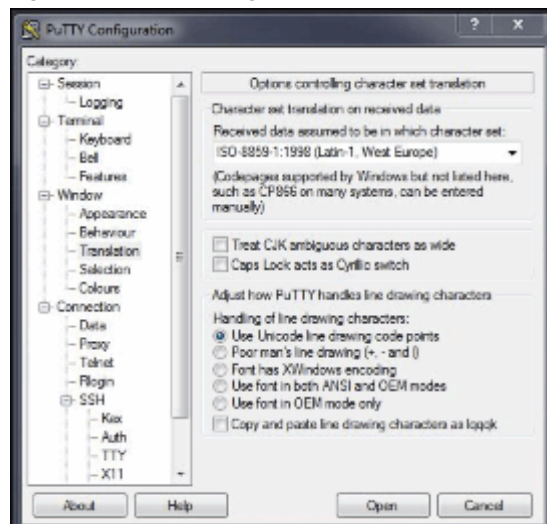

- Under 'Translation' tab.
- Set 'Received data assumed to be in which character set:' to 'ISO-8859-1:1998 (Latin-1, West Europe)'.
- Set 'Handling of line drawing characters:' to 'Use Unicode line drawing code points'.




Fig. 118: PuTTY Configuration - Selection

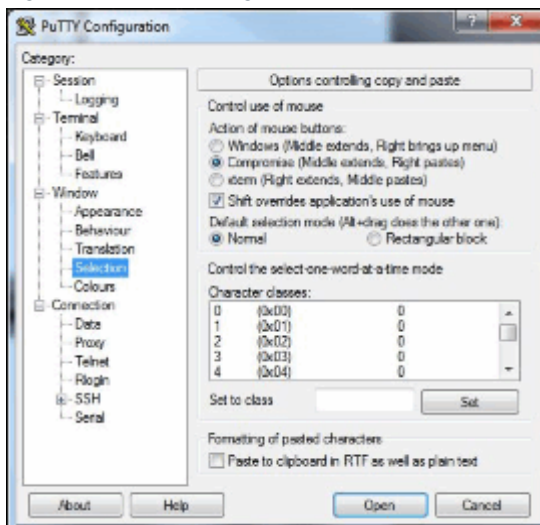


- Under 'Selection' tab.
- Set 'Action of mouse buttons:' to 'Compromise (Middle extends, Right pastes)'.
- Tick 'Shift overrides application's use of mouse'.
- Set 'Default selection mode (Alt + drag does the other one):' to 'Normal'.

Fig. 119: PuTTY Configuration - Colours

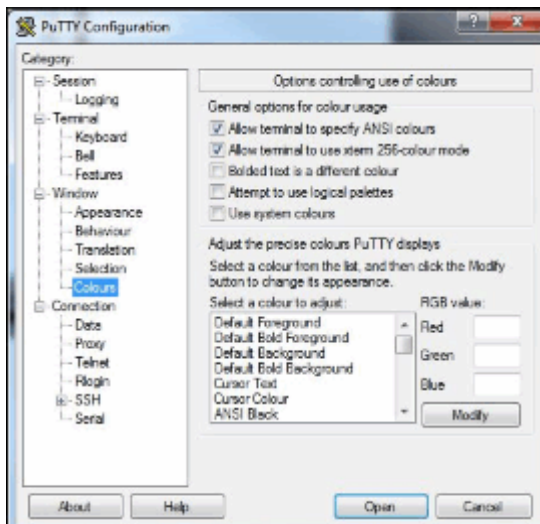


- Under 'Colours' tab.
- Tick 'Allow terminal to specify ANSI colours' and 'Allow terminal to use xterm 256-colour mode'.

Fig. 120: PuTTY Configuration - Connection

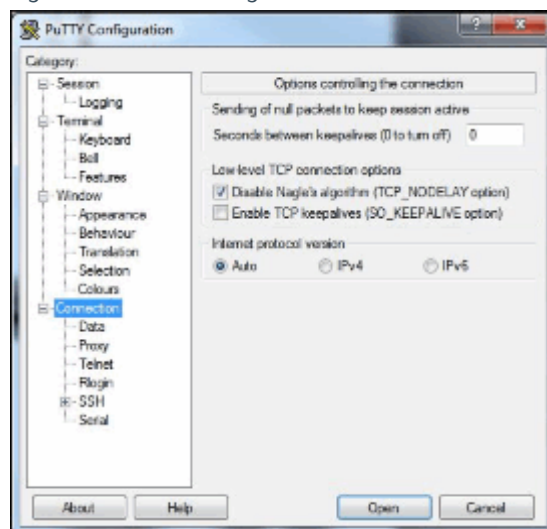


- Under 'Connection' tab.
- Set 'Low-level TCP connection options' to 'Disable Nagle's algorithm (TCP\_NODELAY option)'.
- Set 'Internet protocol version' to 'Auto'.

Fig. 121: PuTTY Configuration - Data

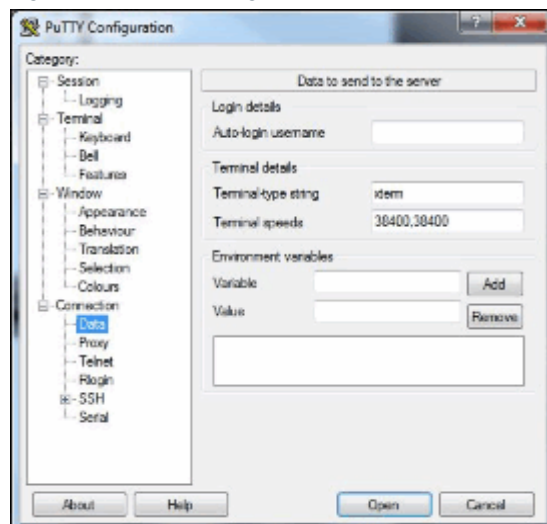


- Under 'Data' tab.
- Under the 'Terminal-type string' ensure it states 'xterm'.
- Under the 'Terminal speeds' ensure it states '38400.38400'.

Fig. 122: PuTTY Configuration - Proxy

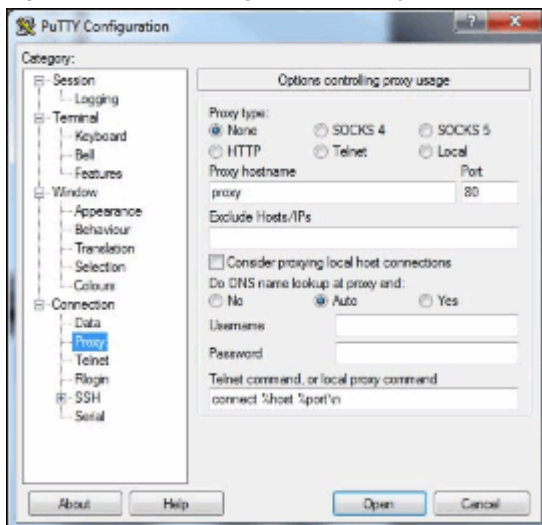


- Under 'Proxy' tab.
- Set 'Proxy type:' to 'None'.
- Under the 'Proxy hostname' ensure it states 'proxy'.
- Under the 'Port' ensure it states '80'.
- Set 'Do DNS name lookup at proxy end:' to 'Auto'.
- Under the 'Telnet command, or local proxy command' ensure it states 'connect %host %port\n'.

Fig. 123: PuTTY Configuration - Telnet

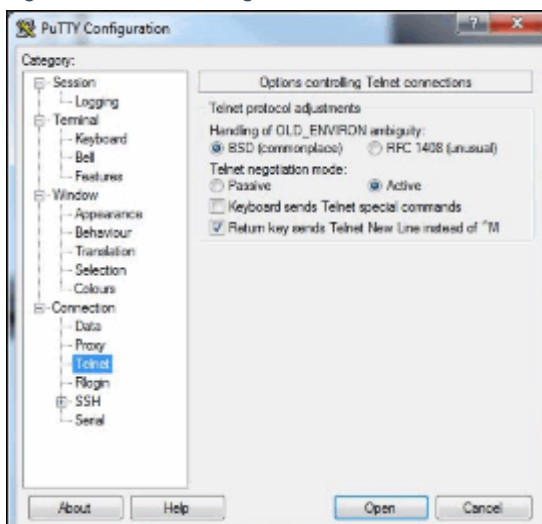


- Under 'Telnet' tab.
- Set 'Handling of OLD\_ENVIRON ambiguity:' to 'BSD'.
- Set 'Telnet negotiation mode:' to 'Active'.
- Tick 'Return key sends Telnet New Line instead of ^M'.

Fig. 124: PuTTY Configuration - Rlogin

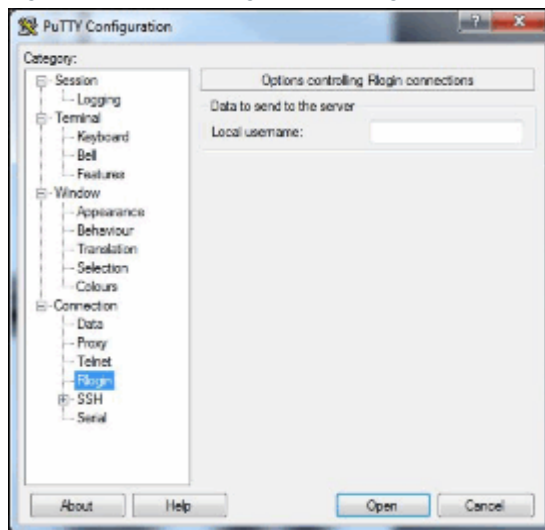


- Under 'Rlogin' tab.
- Nothing is selected.

Fig. 125: PuTTY Configuration - SSH

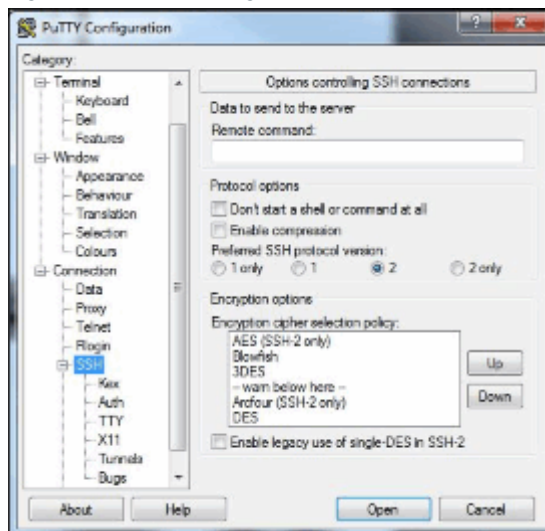


- Under 'SSH' tab.
- Set 'preferred SSH protocol version' to '2'.

Fig. 126: PuTTY Configuration - Kex

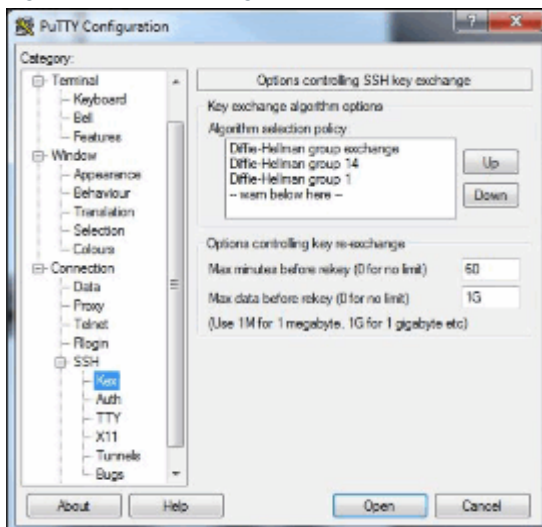


- Under 'Kex' tab.
- Under the 'Max minutes before rekey (0 for no limit)' ensure it states '60'.
- Under the 'Max data before rekey (0 for no limit)' ensure it states '1G'.

Fig. 127: PuTTY Configuration

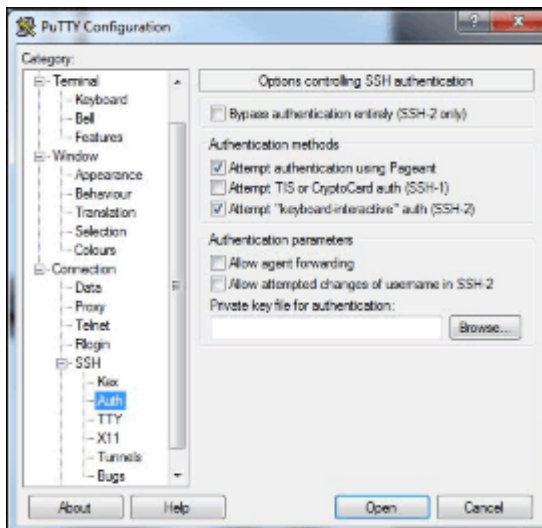


- Under 'Auth' tab.
- Tick 'Attempt authentication using Pageant' and 'Attempt "keyboard-interactive" auth (SSH-2)'.

Fig. 128: PuTTY Configuration - TTY

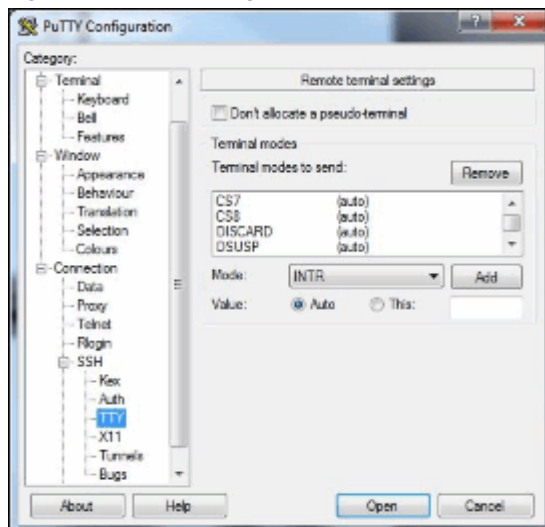

- Under 'TTY' tab.
- Set 'Value' to 'Auto'.


Fig. 129: PuTTY Configuration

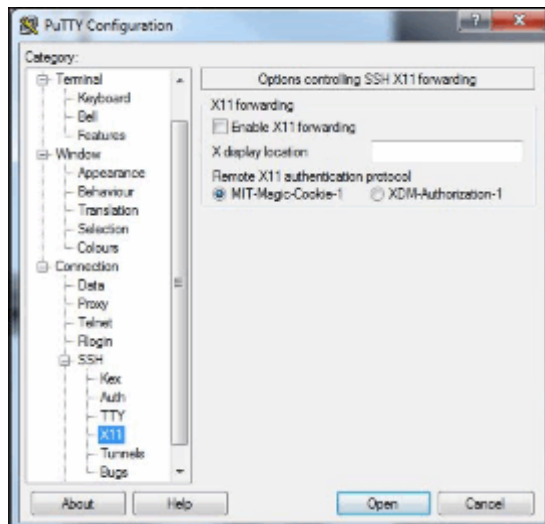

- Under 'X11' tab.
- Set 'Remote X11 authentication protocol' to 'MIT-Magic-Cookie-1'.




Fig. 130: PuTTY Configuration - Tunnels

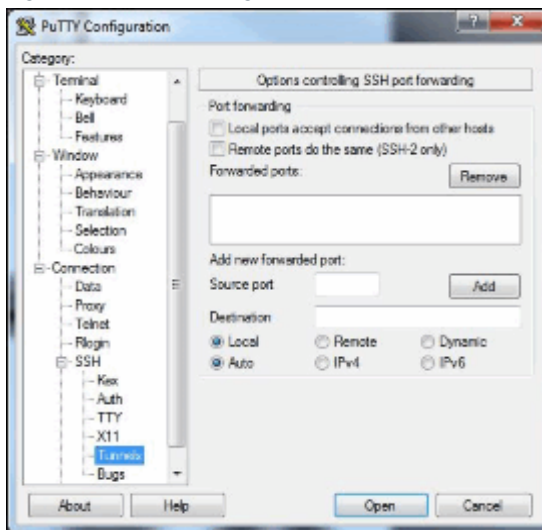


- Under 'Tunnels' tab.
- Set 'Destination' to 'Local' and 'Auto'.

Fig. 131: PuTTY Configuration - Bugs

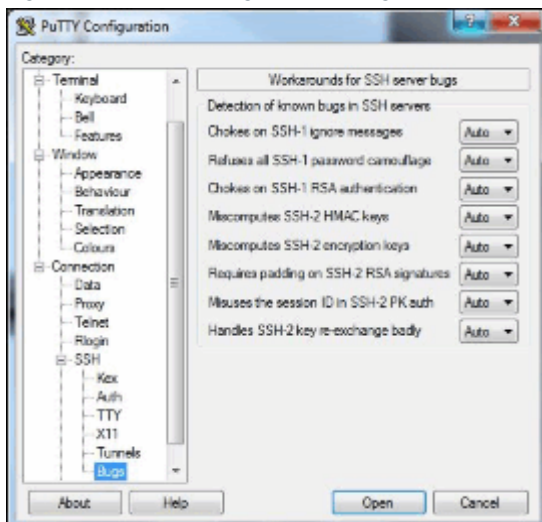


- Under 'Bugs' tab.
- Ensure all settings are set to 'Auto'.

Fig. 132: PuTTY Configuration - Serial

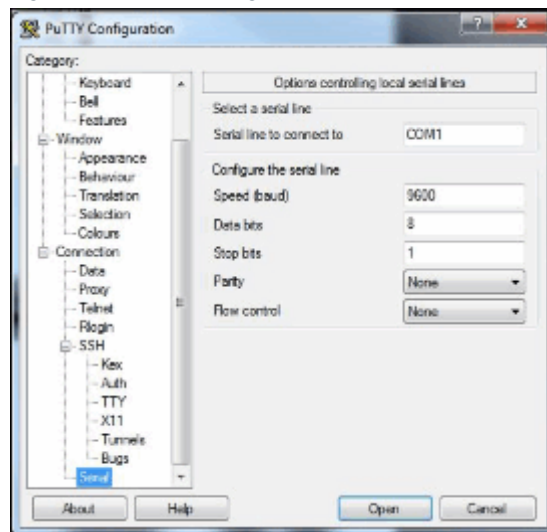

- Under 'Serial' tab.
- Under the 'Serial line to connect to' ensure it states 'COM1'. Only state this if the com port that is being used. Change as required.
- Under the 'Speed (baud)' ensure it states '9600'.
- Under the 'Data bits' ensure it states '8'.
- Under the 'Stop bits' ensure it states '1'.
- Under the 'Parity' ensure it states 'None'.
- Under the 'Flow control' ensure it states 'NONE'.

Fig. 133: PuTTY Configuration - Session

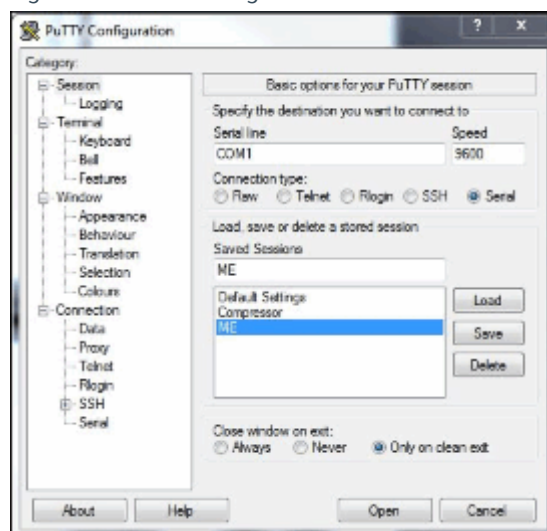

- Under 'Session' tab.
- Under 'Saved Sessions' type in 'ME' or any name which you want to save these settings under.
- Click the 'Save' button.
- Click the 'Open' button.
- PuTTY software will now open and run.

In future requirements to use PuTTY follow this short sequence.

- Open the PuTTY software.
- Under the Session tab select the session required. In the case above the session is called ME.
- Click the load button.
- Click the 'Open' button.
- PuTTY will now open and run.

#### 4.1.4.5 Ramp Cables

- Ramp cables **must not be fitted** until the **Transport test** has been **successfully** completed.

## 4.1.4.6 Initialising the 3600 MPSU



Ensure that there is a valid helium level reading before switching on the 3600 MPSU.


**CAUTION**

Risk of quench.

Failure to observe the following process may cause the magnet to quench.

- » Do not perform any other software sequence or downloads on the Host or Service laptop during ramping.

Fig. 134: Transport test tool

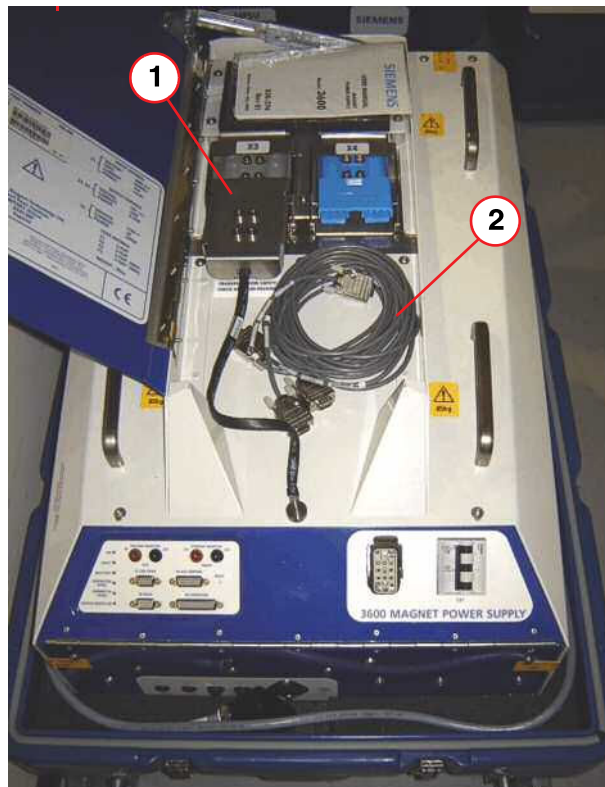

- (1) Transport test plug (fitted)
- (2) RS232 cable

The **transport test tool** consists of 2 x 470 ohm resistors housed in a grey Anderson plug.

The center tap of the resistors is connected to a flying lead which is permanently bolted to the 3600 MPSU chassis.

The MPSU is supplied **with the test plug (X3)** already connected. (Pos 1)

When the power supply arrives at site, the **transport test plug X3** should already be fitted. (→ 1/ Fig. 134 Page 115). If not, connect it into **X3**. The test must be **passed** before it is removed and the **ramp cables** connected.

- Ensure the **RS232 cable** is connected between the **X6 connector** on the 3600 MPSU and the **laptop**.
- Ensure the **Terminal emulation software** is set up and running.



The Transport test is only required for the first time the MPSU is used on a new installation.

- Ensure that **EPC circuit breaker (F53)** is set to the **ON position**.
- **At the MPSU**
  - set circuit breaker **CB1** to the **ON position 'I'**
  - check the **green LED** indicator is '**ON**' and **lit**.
- The '**boot**' screen is displayed.

Fig. 135: Terminal boot-up screen

```
MPS3600 Serial Connection
Magnet Power Supply 3600 551-362 VAOSK 3rd November '11
Enter 3 characters for your initials...
pns
Detecting the magnet serial number....
Magnet serial number received 99000005
Detecting the magnet type.....
Magnet has been identified as OR99
Checking magnet is supported...
Magnet is supported, getting it's data...
Requesting magnet current..
Magnet current received 466.70A
Requesting Helium level.
Helium level OK 59.9%
CCR supported
```



If the screen does not appear, make sure that the correct COM port and settings have been selected. If they have, turn the MPSU off and then back on, using the mains circuit breaker CB1. (→ 1/ Fig. 99 Page 97)



Reboot the scanner with the MPSU switched on.

A Reboot is necessary to update the loadware and to download the ramp logfile to the HOST.

The logfile is saved to C/MedCom/log directory.

Always reboot the scanner after the MPSU is switched off and on again. Rename the logfile, e.g. "MPSRampLog1.dat" after the ramp is finished.

If not renamed, the file will be overwritten with the next ramp logfile.

Repeat booting of the scanner if following warning message still persists: "MPS not operational"

### Mandatory Transport test



During transport, it is possible that the 3600 MPSU may have suffered internal damage due to rough handling.



The mandatory Transport test was introduced from 3600 MPSU software revision level 17. (→ Transport Test Procedure / Page 117)

For units prior to this revision, continue with (→ Ramp lead connection / Page 121).

This test only needs to be completed the first time the 3600 MPSU is used on site.



This test must be performed PRIOR to connecting the ramp leads and using the 3600 MPSU on the magnet.

#### 4.1.4.7 Transport Test Procedure

When prompted, enter your **initials**.


#### If transport test is NOT required

- Turn off the 3600 MPSU via **CB1**.
- **Connect** the **ramp leads** from the **Filter panel** to **X3** and **X4**: (→ Ramp lead connection / Page 121)
- Then refer to **Run to field**.

#### Transport test required

- The 3600 MPSU will now retrieve the **magnet type**, **magnet serial number**, **Last current set**, **He-level**, **date** and **time** from the **MSUP**.
- The following screen will prompt '**Perform Transport Test**'.
  - » Type **Y** and **<ENTER>**.

Fig. 136: Transport Test - prompt

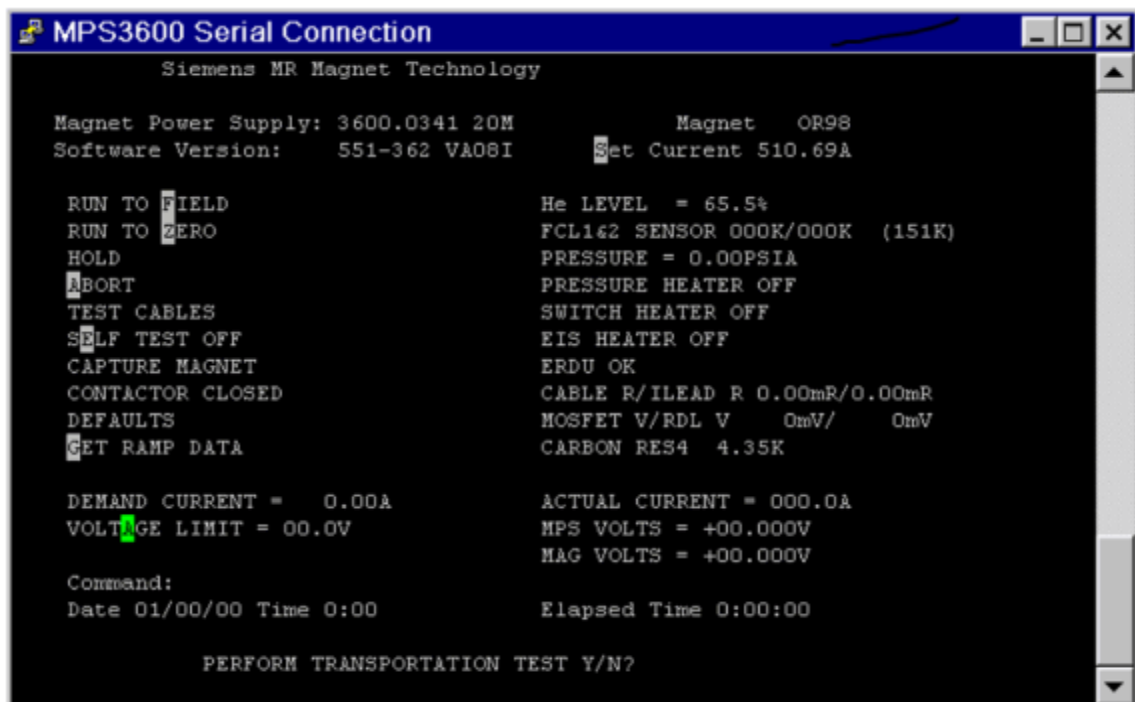

- The 3600 now performs a **transit test**.
- If the test is **successful**, the power supply will perform an **internal self-test**.
- The self test **LED** indicator will **light**.
- Turn off the 3600 MPSU via **CB1**.
- **Connect** the **ramp leads** from the **Filter panel** to **X3** and **X4**: (→ Ramp lead connection / Page 121)
- Then refer to **Run to field**.

#### Failure of the Transport test

**Failure** of the **Transport test** is displayed briefly in the event of the test detecting an **internal fault**. (→ Fig. 137 Page 119).

This is followed by a request to **check the connections** and **repeat** the test. (→ Fig. 138 Page 119)



After 3 successive Transport tests have FAILED, further operation of the 3600 MPSU is prevented. Return the MPSU for repair and order a replacement unit.

Line power to the 3600 MPSU must be switched OFF then ON to enable further tests to be performed.



If all attempts fail, the 3600 MPSU must be returned for repair, and a replacement unit obtained.



Fig. 137: Transport Test failure


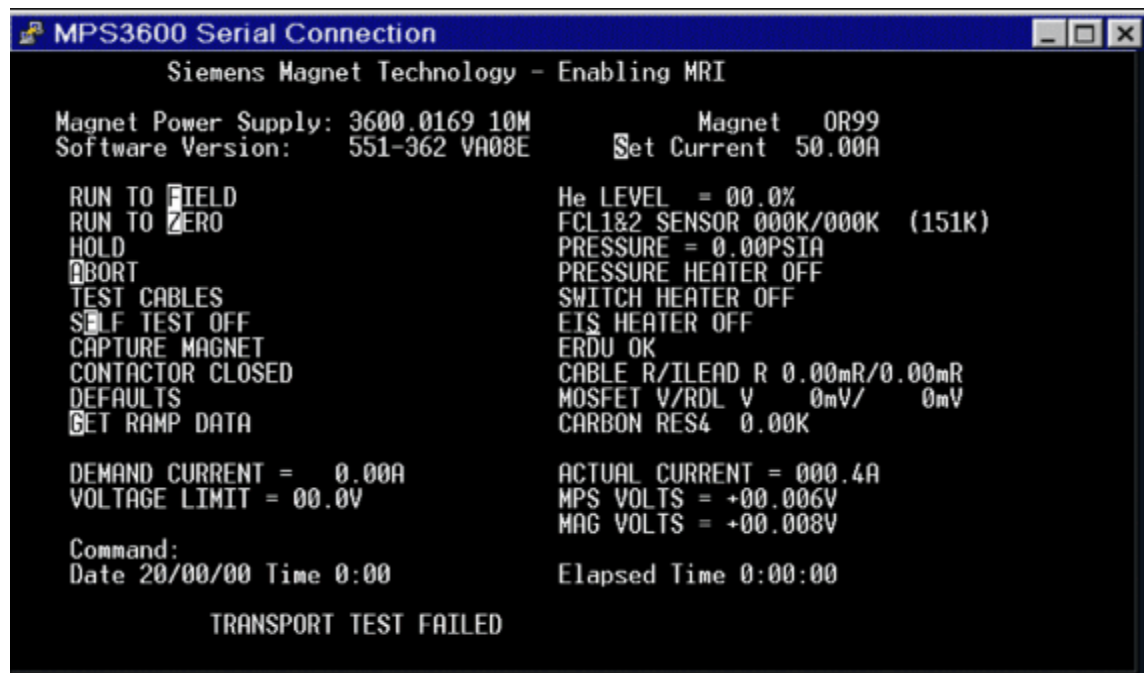


Fig. 138: Transport Test - verify test connector

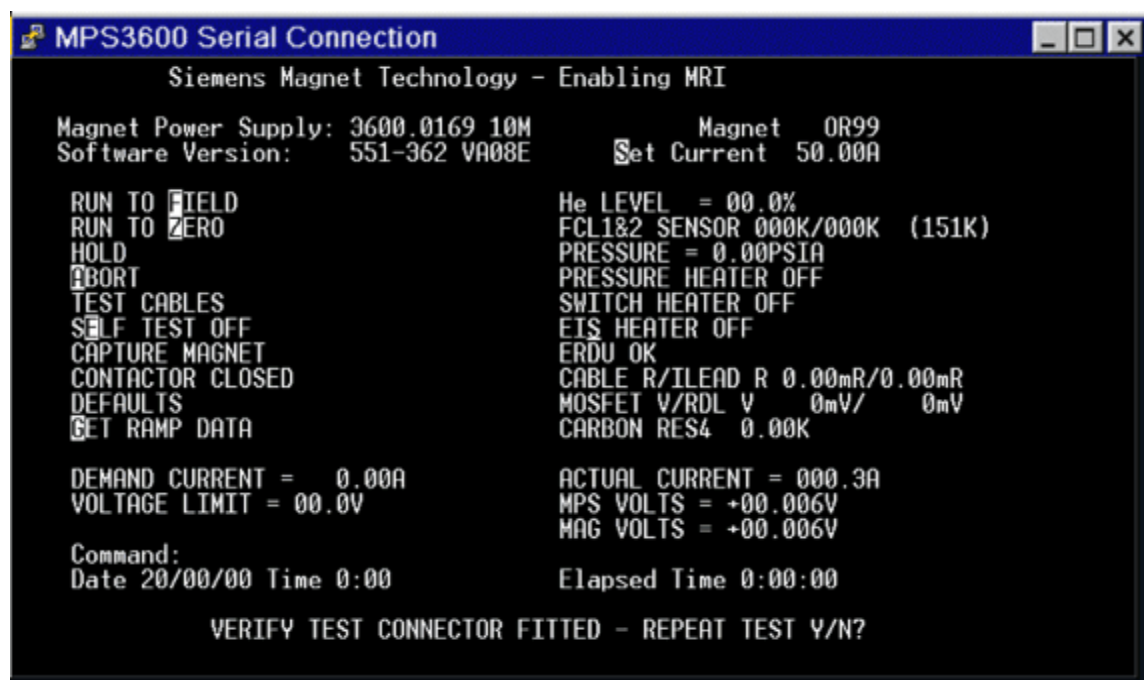


### Internal self-test

During the self-test, the **LED indicator** indicator will light during the test. The self-test runs up to 700A with the output shorting contactor closed.


This enables the voltage across the internal Run-down Load and MOFSET switch to be measured and displayed on the screen.

Fig. 139: Self test

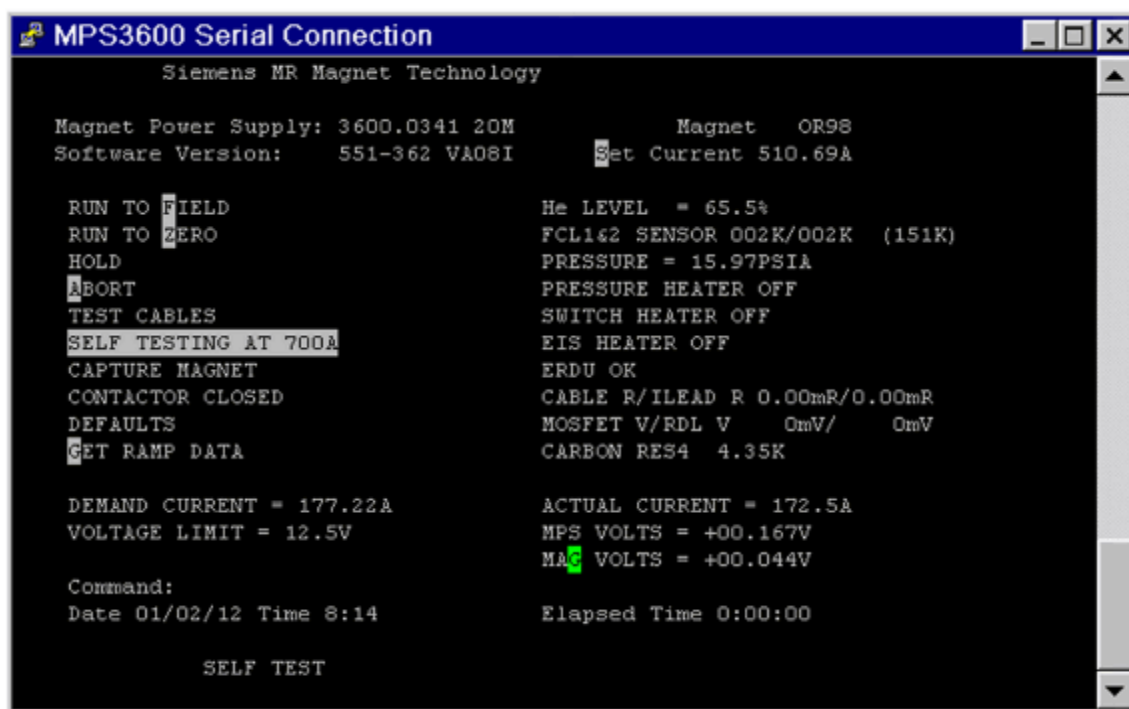

- When the self-test is complete, the **Operational screen** is displayed.

Fig. 140: Operational screen

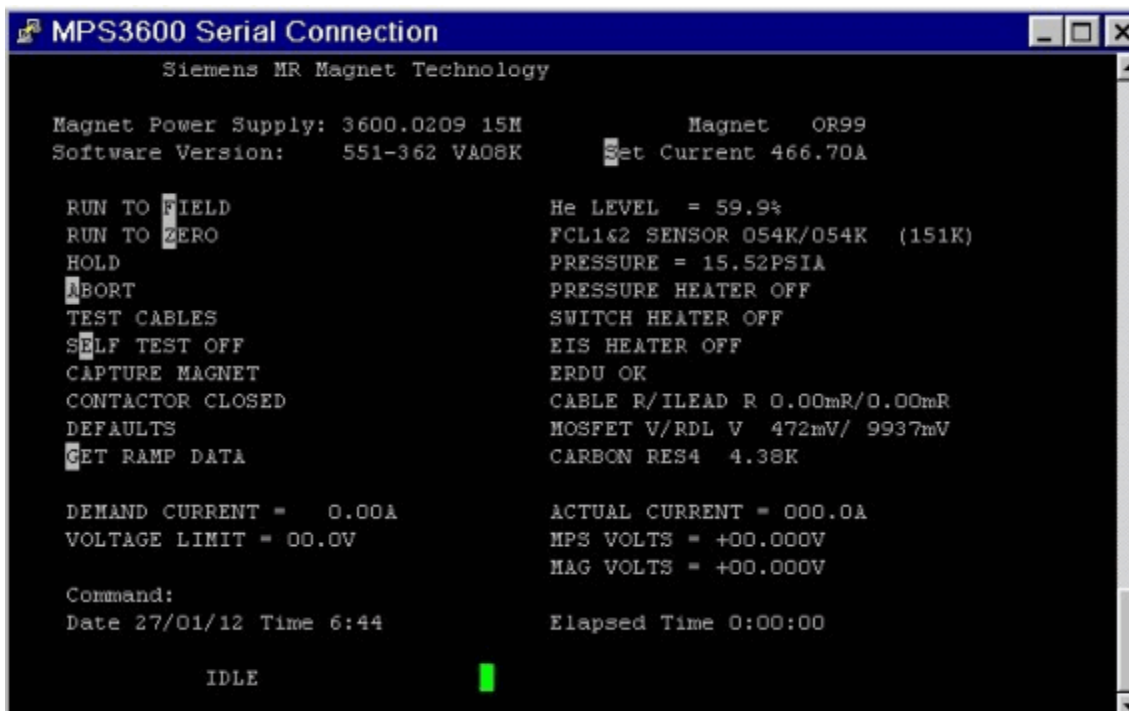

- When packing the MSPU after use, ensure that the **Transport test plug** is refitted, ready for the next site.

#### 4.1.4.8 Ramp lead connection


After the Transport test has been successfully completed, (if applicable), the following connections must be made/changed **prior to ramping the magnet**.



There is no need to switch off the 3600 MPSU at this point.

- **Remove the Transport test plug.**

Fig. 141: Equipment Room - locking plate

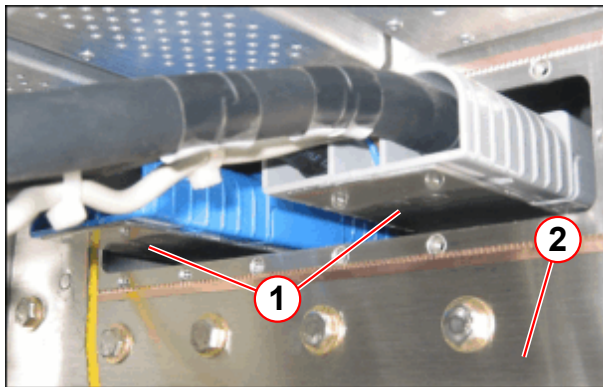

- (1) Locking plate
- (2) RF filter plate

#### Equipment room

- **Connect the 3600 MPSU ramp leads to the magnet leads at the RF filter plate.**


#### Magnet room

- Prevent the **ramp leads** from being accidentally **disconnected**.
  - » **Fit locking bolts and wing nuts.**




Ramp leads are heavy - to prevent stress on the connection, use tie wraps to support the weight if locking bolts are missing.

Fig. 142: Magnet room - locking plate (viewed from below)

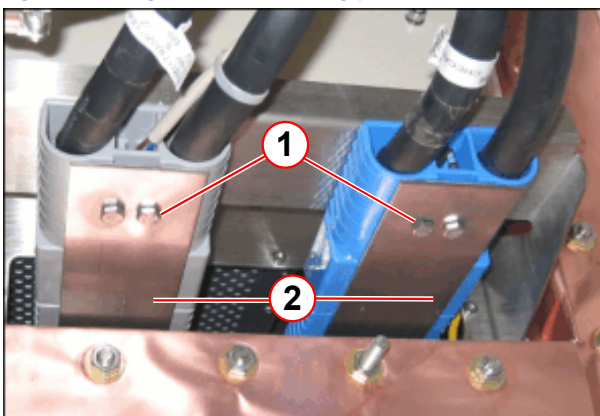

- (1) Locking bolts
- (2) Locking plate

Fig. 143: Magnet room - locking plate (viewed from above)

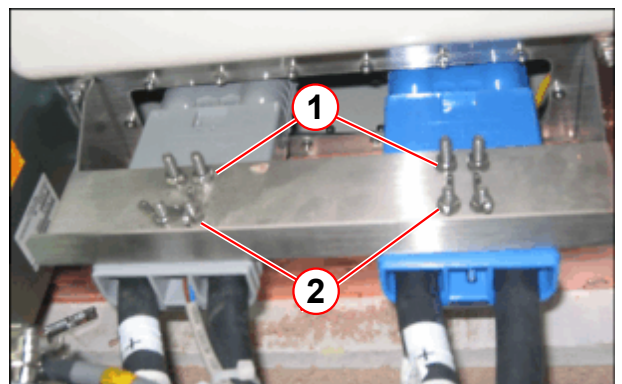

- (1) Anderson fixing bolts
- (2) Wing nuts

Fig. 144: MPSU - ramp cables fitted



- Connect the two **3600 MPSU ramp leads** to Anderson connectors **X3** and **X4** on the top of the **3600 MPSU**.
- The screen will prompt '**ARE RAMP CABLES FITTED**' Y <ENTER>
- The 3600 MPSU now performs an **internal self-test**. (→ Internal self-test / Page 119)

**CAUTION**

Incorrect locking of the ramp leads.

If these instructions are not observed, the magnet and the 3600 MPSU can be damaged.

- » Connect and secure the leads.

### 4.1.5 Run to Field

#### 4.1.5.1 Co-Siting of Magnets

Co-siting describes the ramping requirements where two MR systems with the same operating frequency are sited next to each other in adjacent rooms.



If two systems are installed next to each other in adjacent rooms:

- One magnet must be ramped to the upper frequency limit.
- The other magnet must be ramped to approx. 100 kHz less at the lower frequency limit.



Shimming is only optimized with both magnets ramped up.

Neither magnet is ramping, while the other is running applications. If one magnet is switched off, the other magnet has to be re-shimmed.

#### 4.1.5.2 Run To Field Prerequisites:

- **Wait at least 1 hour** after completing any previous ramp, or aborted ramp (including lead checks).



Following a helium fill, wait at least 2 hours before ramping.

- **Helium** level must be **> 40 %** for the **OR98** and the **OR103**. Or **> 60 %** for the **OR99** for all ramps after the initial ramp.



Initial installation ramp helium level must be **≥ 60 %** for the **OR98** and **≥ 65 %** for the **OR99** and the **OR103**.

- **At the MPSU**
  - set circuit breaker **CB1** to the **ON** position '**I**'
  - check the **green LED** indicator is '**ON**' and **lit**.
- The '**boot**' screen is displayed.
- When prompted, enter your **initials**.

#### 4.1.5.3 Setting the magnet current

The last ramp current will be stored in the MSUP/.iButton and will be displayed on the screen when the 3600 MPSU is switched on. Check the **actual magnet current** required. This is recorded on the **Magnet Acceptance Document**. To change the current:-

- **Enter** the command: **S<space>magnet current<ENTER>**.

#### 4.1.5.4 Run to Field

- **Enter** the command **F <ENTER>** for ramping up the magnet.
  - » **Checks** will be carried out and the 3600 MPSU will then perform an **automatic ramp** to field.

The sequence of screen displays viewed during the ramp up is as follows:

- **Cable test.**
- The MPSU prompts the question '**Has magnet been filled in last 2 hrs Y / N?**'
  - » Enter Y / N **<ENTER>**
- Magnet filled **Yes**
  - » The power supply prompts - **The magnet been filled in the last 2 HRS.**
  - » Answer **Y**.



- Cable test.

Fig. 145: Helium fill check

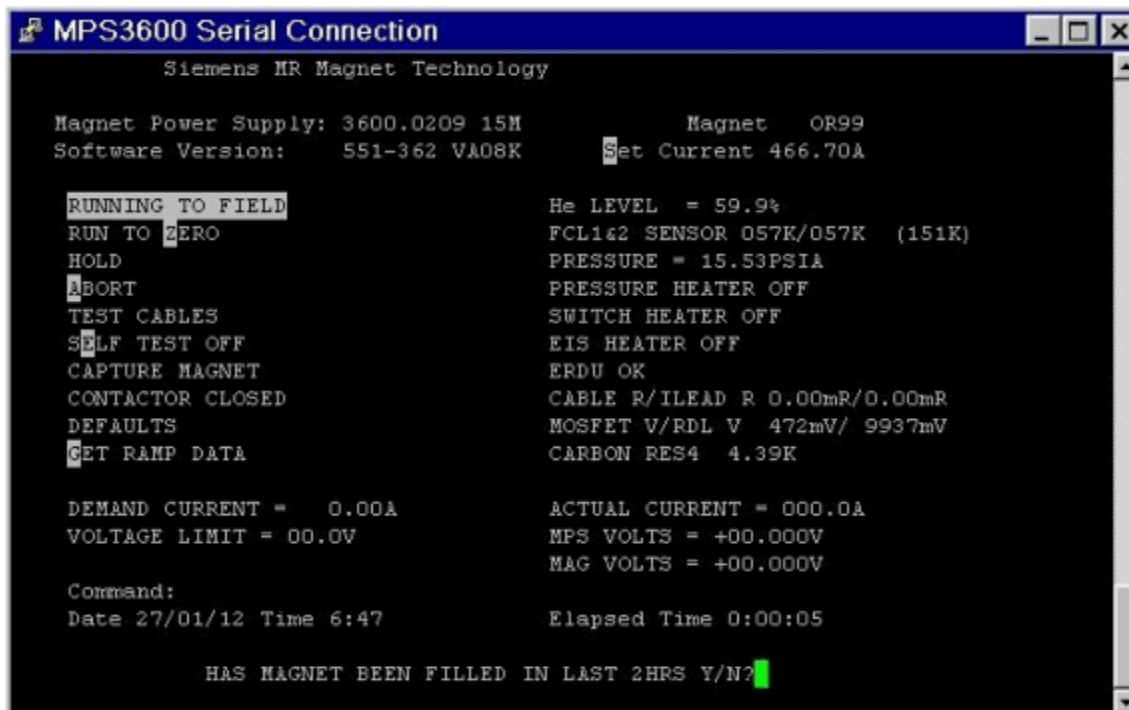


Fig. 146: Pressure heater ON (pressure rising)





- Pressure rising - waiting for **10K drop**.



Pressure raising time will vary, depending on magnet pressures and temperatures.

Maximum 10K drop delay time is:

OR98 = 30 mins


OR99/OR103 = 40 mins



If the 3600 MPSU software times out during the 10K drop the ramp has failed. Abort ramp and wait one hour before restarting.

- Waiting vent delay (5 mins).
- Ramp to field (set current) after performing the overfield algorithm.

Fig. 147: Ramping to field

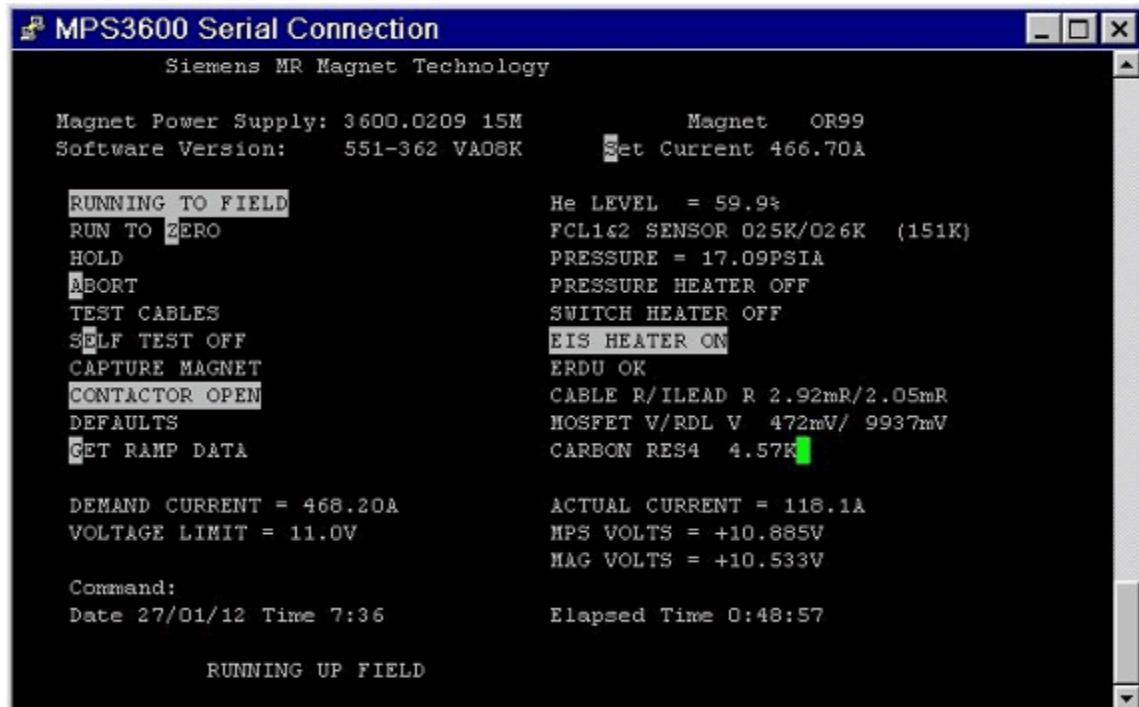

- Field delay (15 seconds).
- Waiting switch delay 2 mins (switch heater off).
- Running lead to zero - EIS delay wait.

- Idle wait - magnet is **at field**.

Fig. 148: Magnet at field (leads down)



- Ensure the **actual current** is equal to zero.
  - » Turn 3600 MPSU **OFF** at **CB1**.



If performing a shim plot:

- Reconnect cable W10220 / X12 / MPSU to W10220 / X12 on top of EPC.
- Restart Software

#### 4.1.6 Run to Zero

##### 4.1.6.1 Run to zero prerequisites:

- Wait at least 1 hour after completing any previous ramp, or aborted ramp (including lead checks).



Following a helium level fill, wait at least 2 hours before ramping.

## 4.1.6.2 Run to Zero

- Ensure the ramp leads are connected to the **3600 MPSU**.




It is not necessary to set the ramp current.

- **At the MPSU**
  - set circuit breaker **CB1** to the **ON** position 'I'
  - check the **green LED** indicator is '**ON**' and **lit**.
- The '**boot**' screen is displayed.
- When prompted, enter your **initials**.
- Ensure the 3600 MPSU is in **Idle mode**.
- Enter the command **Z<ENTER>** for ramping down the magnet.
  - » Checks will be carried out and the 3600 MPSU will then perform an automatic ramp to **zero**.

The sequence during the ramp to zero is as follows:-

- Cable test.

Fig. 149: Helium fill check

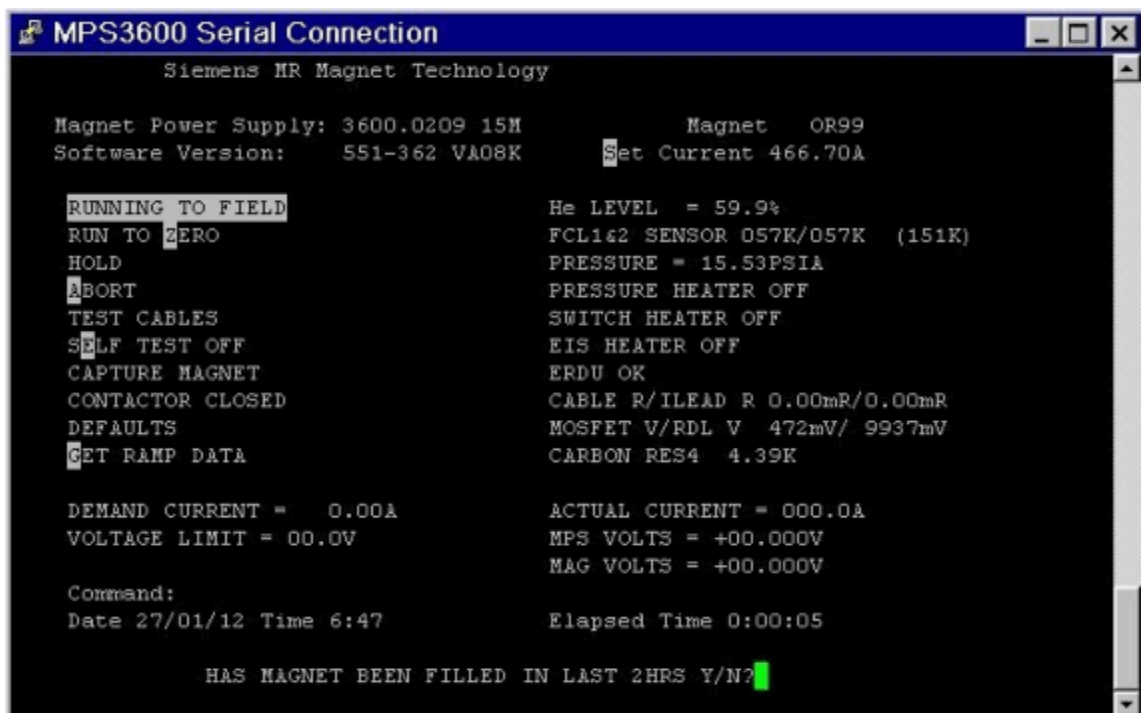

Magnet filled **Yes/No?**


- » The power supply prompts - **Has the magnet been filled in the last 2HRS**
- » Answer **Y / N**

The MPSU unit now controls the run to zero sequence with the following steps:

- 10K drop, 5 minute cooling time.
- Switch heater on signal to MSUP.

- Magnet capture.
- Magnet current running to zero.
- EIS delay wait (after contactor closes).

Fig. 150: Idle wait - Magnet is at Zero

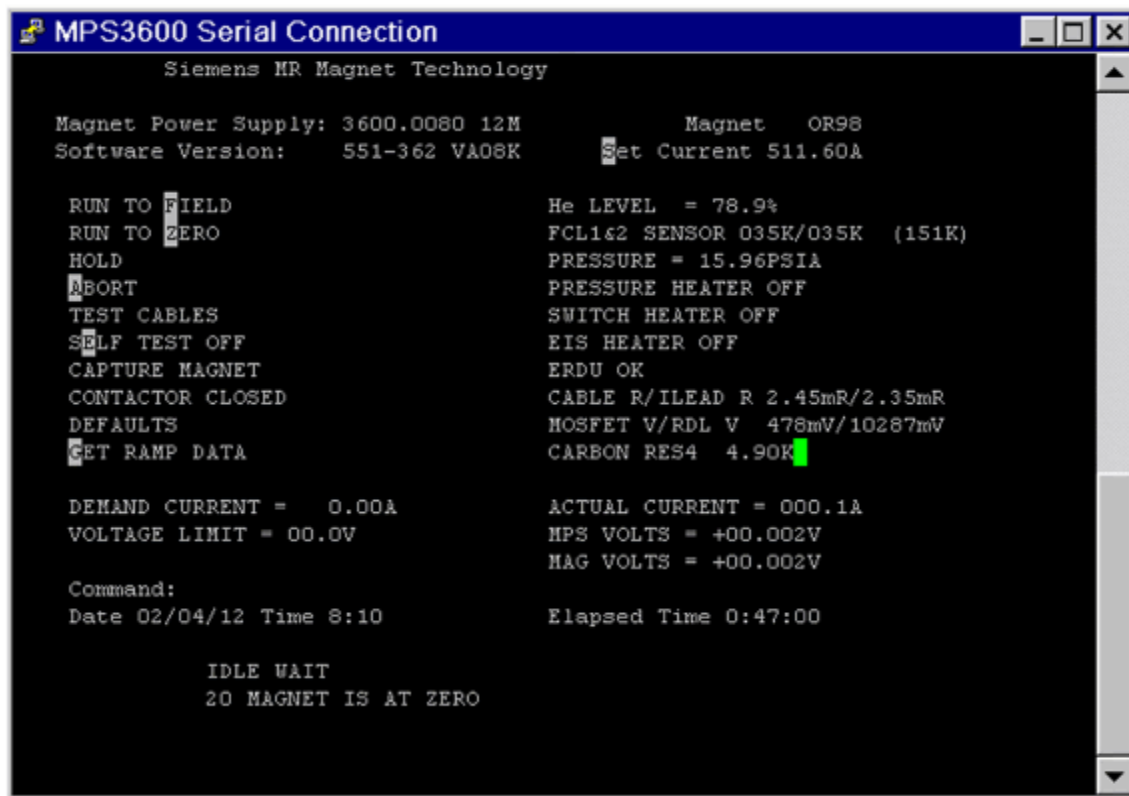

- Idle wait - 20 magnet is **at zero**.
- When the **ramp is complete**,
  - **Download** the ramp logs (→ Downloading the 3600 MPSU log file / Page 130)
  - **Turn off** the 3600 MPSU at circuit breaker switch (**CB1**) - refer to (→ Switching off the MPSU / Page 130).
  - **Disconnect the ramp leads** from the 3600 MPSU.

#### 4.1.7 Ramping the magnet to the required frequency

The magnet must be ramped to the upper limit of the frequency limit. This will require making small changes to the magnetic field current either up or down. The procedure is the same as ramping the magnet.

##### 4.1.7.1 Magnet current setting OR98

Refer to the **Magnet Acceptance Document/MSUP** to obtain the operating **current value**.



OR98 range of frequency:

Maximum frequency: 63,700 000 MHz + 0 kHz

Minimum frequency: 63,500 000 MHz

Use the mathematical formula below to calculate the current to determine the **centre frequency**.

Between ramps, the magnet system's frequency will slowly reduce over a period of time. To increase the period between ramps, aim to achieve for the **upper sector** of the frequency window.

**OR98**

$$I_{\text{final}} = (63,680\,000 \text{ MHz} \times I_{\text{actual}}) / \text{Freq actual}$$

#### 4.1.7.2 Magnet current setting OR99 / OR103

Refer to the **Magnet Acceptance Document/MSUP** to obtain the operating **current value**.



OR99 / OR103 range of frequency:

Maximum frequency: 123,265 000 MHz + 0 kHz

Minimum frequency: 123,100 000 MHz

Use the mathematical formula below to calculate the current to determine the **centre frequency**.

Between ramps, the magnet system's frequency will slowly reduce over a period of time. To increase the period between ramps, aim to achieve for the **upper sector** of the frequency window.

**OR99 / OR103**

$$I_{\text{final}} = (123,260\,000 \text{ MHz} \times I_{\text{actual}}) / \text{Freq actual}$$

#### Magnetom Skyra with Multinuclear Spectroscopy Option (MNO)



OR99 / OR103 range of frequency:

Center frequency: 123,200 000 MHz

Maximum frequency: 123,265 000 MHz +0 kHz

Minimum frequency: 123,125 000 MHz -0 kHz

Use the mathematical formula below to calculate the current to determine the **centre frequency**.

Between ramps, the magnet system's frequency will slowly reduce over a period of time. To increase the period between ramps, aim to achieve for the **upper sector** of the frequency window.

OR99 / OR103

 $I_{\text{final}} = (123,200\,000 \text{ MHz} \times I_{\text{actual}}) / \text{Freq actual}$ 

### 4.1.8 Work after ramping the magnet

#### 4.1.8.1 Downloading the 3600 MPSU log file

The last ramp data is stored in a log file within the 3600 MPSU, even if the 3600 MPSU is switched off. Data is overwritten when the next ramp is started. It is recommended that after each ramp, the log file is downloaded and stored for future analysis.

The Magnet Peak Current is recorded. If the magnet quenches, the recording will tell you at what current level the magnet quenched.



The ramp log files may be saved locally if the host is operating. However it is still recommended to recover these logs as described manually to ensure data consistency.

To download the log file:

- **G** <Enter>, this will ensure that you are in the **error screen**.
- From the Terminal menu, select:-
  - Transfer.
  - Capture Text.
  - Browse.
  - Select the target directory and enter a filename e.g. **C:\Medcom\log\MR25123\_51099log.txt**.
  - Start.
- **L** <Enter>, the data screen appears.
- A scrolling data screen will appear. Wait until the data scrolling has stopped.  
The data screen contains seven columns. Each line displays the data stored at a specific period of time during magnet ramping.
- From the Terminal menu, select:
  - Transfer.
  - Capture Text.
  - Stop.

#### 4.1.8.2 Switching off the MPSU



Verify that the yellow "CURRENT IN LEADS" LED is OFF when the 3600 MPSU is switched on. (→ 5/ Fig. 102 Page 99)





NEVER attempt to remove the magnet ramp leads when the “CURRENT IN LEADS” LED is flashing.

- Set the mains circuit breaker **CB1** to the **OFF** position. (→ 1/ Fig. 99 Page 97)
- Set the circuit breaker **F53** in the **EPC** to the **OFF** position. (→ 1/ Fig. 105 Page 100)
- **Disconnect** the 3600 MPSU ramp leads from the magnet leads at the RF filter panel.
- The **magnet ramp leads** remain fitted and stored on the lead tray in the magnet room.
- **Re-fit** the equipment room **adapter plate** using the securing screws.



Connect the Anderson connectors (blue/blue and grey/grey).

Fig. 151: RF filter plate



Fig. 152: Adapter plate (view from Equipment Room)



- Reconnect cable W10220 / X12 / MPSU to W10220 / X12 on top of EPC.
- Restart the Software.

## 4.1.9 Final work steps

### 4.1.9.1 Save the magnet current to the HOST

- Access the service menu in the **SeSo at MRC (MRAWP) Host**.
- Select menu item **“Magnet & Cooling”**.
- Select **“Magnet Setup”** in the **“Action”** menu.
- Click **“Download”** to save the magnet current value on the **HOST**.

Fig. 153: Magnet setup

Siemens Med Service Software - Windows Internet Explorer

Magnet & Cooling | Terminal | Event Log | Reports | Help

**HOME**

**Magnet & Cooling**

- Actions
- Cooling Status
- Magnet Status
- Compressor Status
- Cooling Setup
- Magnet Setup**

### Magnet & Cooling - Magnet Setup

**Magnet**

|                                     |                      |          |         |
|-------------------------------------|----------------------|----------|---------|
| Helium Level Probe 1                | Min. (0%)            | 85.1     | Ohm     |
|                                     | Max. (100%)          | 0.2      | Ohm     |
| Helium Level Probe 2                | Min. (0%)            | 85.1     | Ohm     |
|                                     | Max. (100%)          | 0.2      | Ohm     |
| Helium Level Warning                | Min.                 | 55       | %       |
| Pressure Heater Resistance          |                      | 20.5     | Ohm     |
| Shield Temperature Warning 50K Link | Max.                 | 74.0     | K       |
| Shield Temperature Warning 50K Bore | Max.                 | 84.0     | K       |
| CCR Calibration                     | Room Temperature     | 0.0      | Ohm     |
|                                     | Nitrogen Temperature | 0.0      | Ohm     |
|                                     | Helium Temperature   | 0.0      | Ohm     |
| Time Setting                        | Daily at             | 02:00    | O'clock |
| MSUP                                | Serial Number        | 1054     |         |
|                                     | Revision Number      | 0        |         |
| Magnet                              | Serial Number        | 21176154 |         |
|                                     | Revision Number      | 0        |         |
| Magnet Current                      |                      | 466.64   | A       |

**MSUP Data**

Downloads MSUP data from MSUP to host [Download](#)

Uploads MSUP data from host to MSUP [Upload](#)

**MSUP Logfiles**

Downloads all MSUP logfiles to the host [Download](#)

☐ Save values as default

Ready: < | Save | Default | Update

#### 4.1.9.2 Warm Turret & Venting

Fig. 154: Warm the turret and venting with a heat gun.



- Ensure that the **turret and venting** remains leak tight and the valves continue to function correctly.
- **Warm** up the **turret and venting** with a heat gun.


#### 4.1.9.3 Pack the 3600 MPSU

Fig. 155: MPSU - ramp cables fitted



- Disconnect the two **3600 extension ramp leads** from the Anderson connectors **X3** and **X4**.
- Re-fit the **Transport Test Plug** into **X3**.
- Disconnect the remaining leads and pack inside the transport case.
- Close the lid and secure by pulling the two silver knobs and turning them through a 1/4 turn.

Fig. 156: MPSU

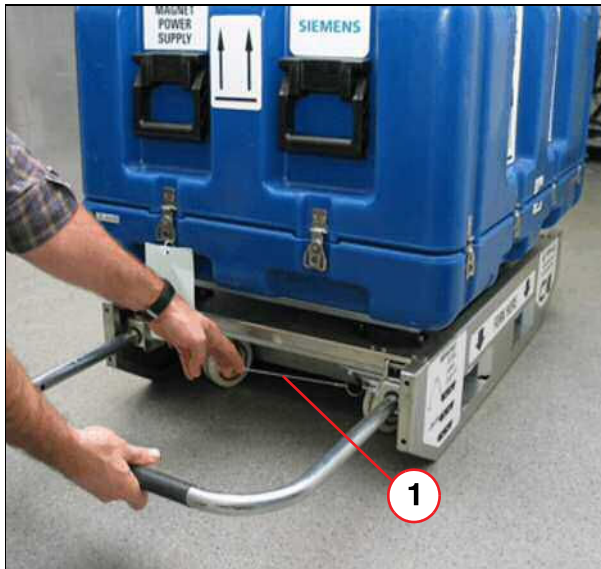

(1) Wire locking strap

- Replace the shipping case lid and secure using the 12 spring clips.
- Pull the wire strap to unlock the handle, and slide the handle out of the recess under the trolley.
- Release the wire strap to lock the handle in the extracted position.
- Raise the handle to lower the wheels and lift the trolley off of the two skids.
- The 3600 MPSU and trolley can now be safely moved.

#### 4.1.10 Mains failure during ramping



#### CAUTION

If the supply to the 3600 MPSU fails while the magnet is ramping through the ramp leads up or down, the magnet current will continue to flow.

Failure to observe the following safety warning may result in physical injury and equipment damage.

- » **Never remove the ramp leads from the magnet while current is in the magnet ramp leads.**

Current is in the ramp cable whenever:

- » The yellow "CURRENT IN LEADS" LED is flashing.



Failure of the 3600 MPSU may prevent this LED from lighting up. So, in the event of a failure (e.g. mains failure), it is essential to measure the voltage across the terminals on the turret.

Any voltage on these terminals other than zero is an indication that there is current in the ramp leads.

- » Magnet was ramping up or down and there was a mains failure.
- » Magnet was ramping up or down and the 3600 MPSU detected a fault such as over temperature, which would shut down the 3600 MPSU.
- » Normal ramp up or down of the magnet.

**CAUTION**

If the magnet is suspected of ramping down through the ramp leads, owing to one of the reasons above.

Failure to observe the following safety warning may result in physical injury and equipment damage.

- » Do not remove the ramp leads. The magnet must be run to zero or persist at field before removing the ramp leads.



If the cause of the shutdown was failure of the mains, nothing can be done until both the 3600 MPSU and the Magnet Supervisory (MSUP) are working again.

When the 3600 MPSU and the MSUP are operational, ramp the magnet to zero.


Do not delay. Ramp down as soon as mains power returns.



Do not wait one hour to ramp down or capture.

#### 4.1.11 iButton function

Fig. 157: i-Button location

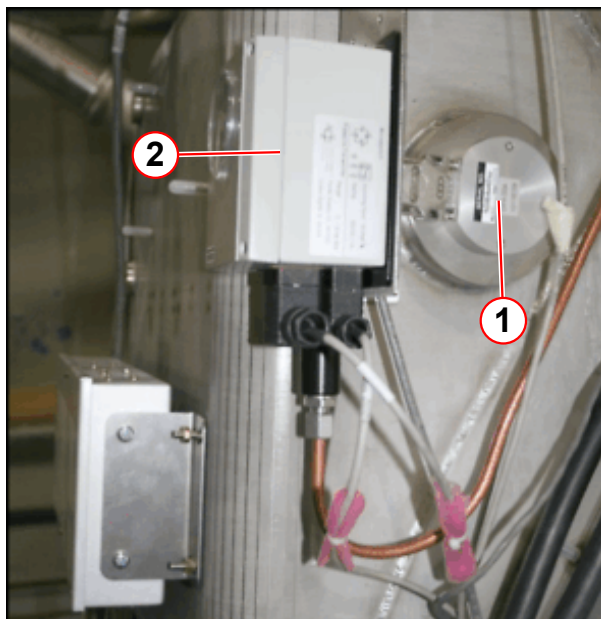

- (1) iButton
- (2) Pressure gauge

The **iButton** is a non-volatile memory device which stores vital magnet data. It is connected to the **MSUP** at **X2** (CH/DIAG/PRESS/ID/VALVE).



## 4.2 What to do after a quench

### 4.2.1 General

If the magnet quenches, a large volume of helium gas is generated within the helium vessel. The quench valve opens, and stay open until the flow of gas subsides. The quench valve re-seals, preventing air (caused by cryo-pumping) from being drawn back into the magnet turret and helium vessel.

The MSUP automatically maintains the system pressure above atmospheric.



Magnet HSC must be notified of a quench (→ Siemens Healthineers Magnet WS CS HSC, <mailto:magnet.wscshsc.team@siemens-healthineers.com>)

Fill out the **Quench Report** - (→ MR Product Information, <https://productinformation.healthcare.siemens.com/ProductInformation.aspx?uuid=7abf6bb6-5daf-440c-9920-fab2d4640084>) - Links - MR - Quality Database (Quench Report )



The magnet helium level can be recovered if the temperature of the CCR is still colder than nitrogen temperature.

The CCR (CARBON RES 1) reading must be in the range of between 1200 and 3400 ohms. Compare to the values on the Magnet Acceptance Document.

Check the helium level at least 1 hour after a quench. If there is 10% helium in the magnet re-fill using standard fill procedures.

To re-fill the OR99 or OR103 magnet from zero helium, or warm ceramic carbon resistor (CCR) temperatures between LN2 and He, perform a BOTTOM FILL procedure.

OR98 magnet requires approximately 1000 to 1500 liters of liquid helium.

OR99 and OR103 magnet requires approximately 2000 liters of liquid helium.

### NOTICE

If the magnet temperature warms up above 77K, a full nitrogen pre-cool is required.

Helium filling will not cool the magnet sufficiently to collect liquid helium.

» A full nitrogen pre-cool will take 8 –10 days.

#### 4.2.1.1 Safety



Post quench procedures must be performed by CSE's that are trained to the MR2MAG-NET level. (or System level trained).





Switch off and isolate the magnet power supply from the mains before further action.



### WARNING

Risk of asphyxiation.

Failure to observe the following safety warning may result in dizziness and loss of consciousness in personnel working with cryogenics.

- » Do not vent helium gas into the magnet room during helium transfer. As an exception, a small quantity of helium gas can be vented to cool the syphon.
- » Ensure the room is ventilated.
- » Ensure that the turret venting kit is connected to the quench pipe.
- » O<sub>2</sub> monitors must be used in the imaging suite.



### CAUTION

**RISK OF ICING!!**

Air ingress will cause INTERNAL ICING OF THE TURRET.

- » Replacement of burst disk(s) must be performed by a qualified CSE.
- » Replace broken or damaged burst disk and O-ring seals IMMEDIATELY.
- » Open the bypass valve to check that there is a flow of gas. If there is no gas flow, there may be an ice blockage. Close the bypass valve.
- » If ice is suspected contact Magnet HSC immediately.



### CAUTION

Cold helium gas.

Personal injury may occur.

- » Protective clothing (long sleeves), cryogenic gloves and eye protection must be worn.


## 4.2.2 Post Quench Steps

Follow the work steps in the table below for the different types of quench.

|        | Type of quench                                                                                                                                                                                                                                                                                                                                                                      |                                                                 |                                                                                                                        |
|--------|-------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|-----------------------------------------------------------------|------------------------------------------------------------------------------------------------------------------------|
|        | Quench during ramping                                                                                                                                                                                                                                                                                                                                                               | Quench during He filling                                        | Persistent mode quench or ERDU button pressed                                                                          |
| Step 1 | <ul style="list-style-type: none"> <li>Isolate the MPSU from the magnet.</li> <li>Download the 3600 MPSU ramp log files.</li> </ul>                                                                                                                                                                                                                                                 | <ul style="list-style-type: none"> <li>Stop filling.</li> </ul> | <ul style="list-style-type: none"> <li>Check that no ERDU buttons are pressed.</li> <li>Re-set ERDU button.</li> </ul> |
| Step 2 | <ul style="list-style-type: none"> <li>Check the liquid helium level.</li> <li>Check the ceramic carbon resistor temperatures at the Host MRC- Magnet &amp; Cooling. (→ Check the Ceramic Carbon Resistor (CCR) Values / Page 142)</li> </ul>                                                                                                                                       |                                                                 |                                                                                                                        |
| Step 3 | <ul style="list-style-type: none"> <li>Order <b>1000 - 2000</b> liters of liquid helium.<sup>1</sup></li> </ul>                                                                                                                                                                                                                                                                     |                                                                 |                                                                                                                        |
| Step 4 | <ul style="list-style-type: none"> <li>Warm up and leak check the quench valve. (→ Checking the Quench Valve and Burst Disk / Page 150)</li> <li>Warm up and leak check the auxiliary vent burst disk. (→ Checking the Auxiliary Burst Disk / Page 154)</li> </ul>                                                                                                                  |                                                                 |                                                                                                                        |
| Step 5 | After 1 hour : <ul style="list-style-type: none"> <li>Download the Magnet Supervisory log files. (→ Log File Capture / Page 141)</li> <li>Helium fill. Fill to recommended helium ramp level.</li> </ul>                                                                                                                                                                            |                                                                 |                                                                                                                        |
| Step 6 | <ul style="list-style-type: none"> <li>Leak check the quench valve and venting.</li> </ul>                                                                                                                                                                                                                                                                                          |                                                                 |                                                                                                                        |
| Step 7 | <ul style="list-style-type: none"> <li>Fill out the Quench report.</li> </ul>                                                                                                                                                                                                                                                                                                       |                                                                 |                                                                                                                        |
| Step 8 | <ul style="list-style-type: none"> <li>Email the (→ Magnet HSC, <a href="mailto:magnet.wscshsc.team@siemens-healthineers.com">mailto:magnet.wscshsc.team@siemens-healthineers.com</a>) <b>with:</b> <ul style="list-style-type: none"> <li>» 3600 Power Supply logs.</li> <li>» Magnet Supervisory logs and screen capture.</li> <li>» Magnet serial number.</li> </ul> </li> </ul> |                                                                 |                                                                                                                        |
| Step 9 | <ul style="list-style-type: none"> <li>Ramp the magnet.</li> </ul>                                                                                                                                                                                                                                                                                                                  |                                                                 |                                                                                                                        |
| Step10 | <ul style="list-style-type: none"> <li>Perform the shim check.               <ul style="list-style-type: none"> <li>» If the shim check exceeds the specifications, a re-shim is necessary.</li> </ul> </li> </ul>                                                                                                                                                                  |                                                                 |                                                                                                                        |

1. The quantity of the order depends on the helium level remaining in the system after a quench, and volume already on site.

Fig. 158: Checking helium level - Post Quench

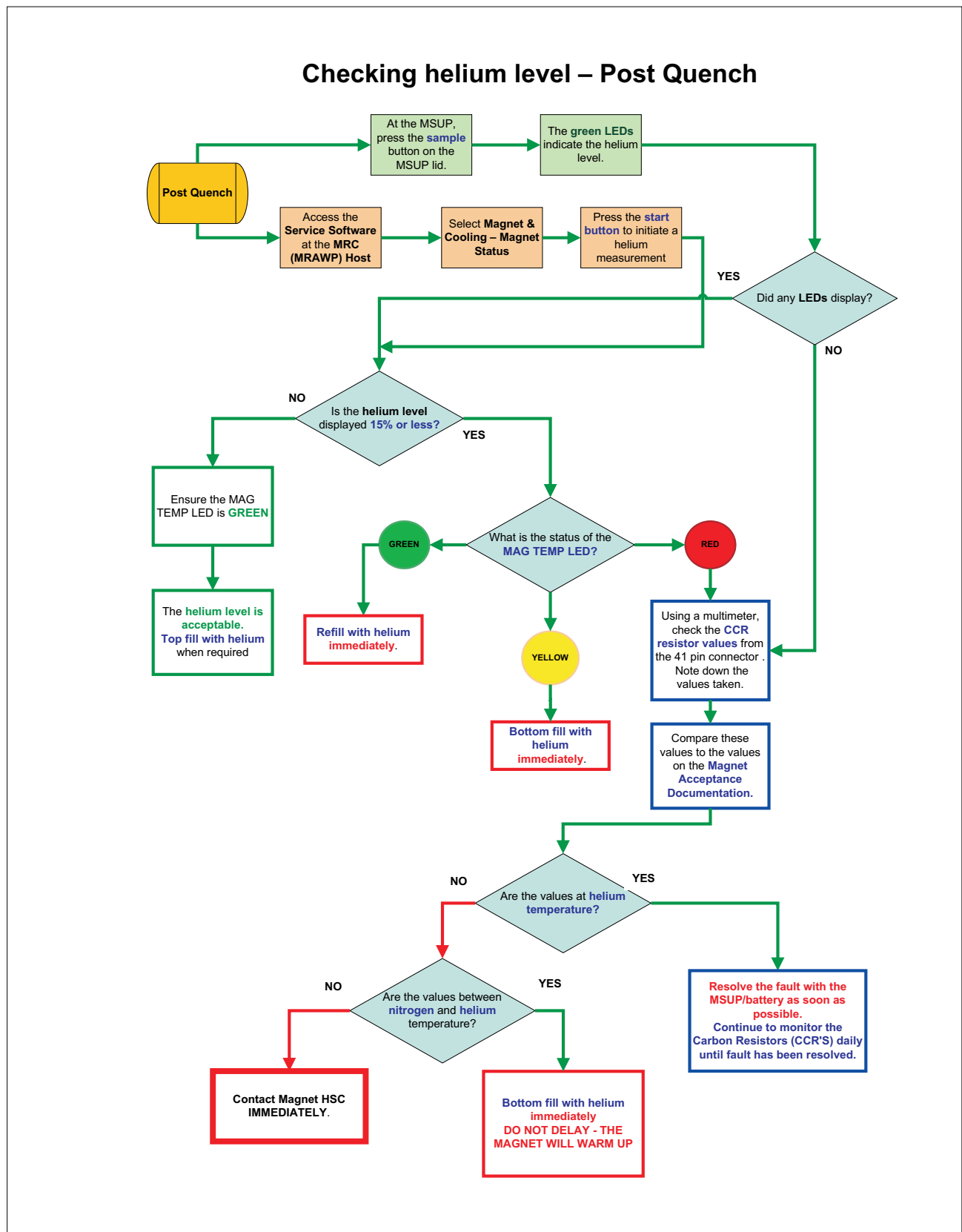

### 4.2.3 Pre-requisites

#### 4.2.3.1 Tools Required

The post quench tools required:

- Heat gun and paper toweling.
- Service tool kit.
- Laptop and RS232 cable.
- O<sub>2</sub> monitor.
- Personal protection equipment. To include cryogenic gloves, safety shoes, glasses and a long sleeved protective clothing.
- Helium leak detector.

Fig. 159: Edwards Leak detector (example)



- Use the Edwards hand-held leak detector (material number **10614850**). This is the recommended equipment for leaks greater than  $1 \times 10^{-4}$

**NOTE:** The hand-held device is calibrated in ml/sec:-

$1 \times 10^{-4}$  mbar l/s =  $1 \times 10^{-4}$  ml/sec at normal pressure and temperature.

#### 4.2.3.2 Spare Parts

Helium leaks from the magnet pressure vessel are generally repaired by replacing a damaged seal.

- A number of seals are supplied in the spares kit. These must be replaced when used on the magnet.
- Venting components will require ordering from CSML.

##### Spare parts for the Quench valve

- Burst disk 15 psi (material number **10128479**).
- Burst disk silicon O-ring, 104.5 ID (material number **10122816**).
- Quench valve high load seal (material number **10128278**).
- Quench valve gasket (material number **10127828**).

##### Spare parts for the Auxiliary Burst disk

- Auxiliary 1.38 bar (20 psi (G)) burst disk (material number **10128480**).
- Auxiliary vent high load seal (material number **10128244**).

- Auxiliary blanks (From supplied spares kit).
  - » Auxiliary vent fiber gasket seal (material number **10128492**).
  - » Red tooling blanks or gasket blanking seal.



Always replace the burst disks and seals used during replacement procedures. They are critical components for reducing the risk of air ingress. A replacement must always be available for immediate use.

#### 4.2.3.3 Unattended Quench

Unattended quench indicators are:

- If the Magnet STOP (ERDU) has been activated:
    - » the yellow WARNING LED lights on the alarm box.
    - » the alarm buzzer is activated.
  - If the helium level is low:
    - » the yellow WARNING LED lights on the alarm box.
    - » the alarm buzzer is activated.
  - If the magnet pressure is very low, the
    - » cold head switches off.
    - » the yellow WARNING LED lights on the alarm box.
    - » the alarm buzzer is activated.
  - Imaging is not possible if there is no magnetic field.
- Press the AUDIO ALARM OFF button on the alarm box to silence any audible alarms.

#### 4.2.3.4 Log File Capture

##### Magnet Supervisory

- The **quench data is lost**, if:
  - a) The system is re-booted within 1 hour after a quench
  - b) The system is turned off within 1 hour after a quench
  - c) The data is downloaded within 1 hour of after a quench
- Download the MSUP log files. (**TSG Downloading MSUP 2306 data to the MRC**).

**If the magnet quenches during ramping:-**

- Download the Magnet Power Supply log file. (**Startup - Downloading MPS Log files**).



The MPSU log files **MUST** be downloaded **BEFORE** the next ramp is performed, otherwise previous data will be overwritten!

- Switch off and disconnect the MPSU from the magnet.

If the magnet quenches during helium fill:

- Stop the transfer.

i

Do not remove or disconnect the syphon as the magnet will be at high pressure.

i

For all types of quench, allow up to 1 hour for the quench pressure to subside.

4.2.3.5 Check the Ceramic Carbon Resistor (CCR) Values

i

Compare the temperatures with the values on the Magnet Acceptance Document.

- Access the Service Software at the MRC (MRAWP) Host.
- Select the menu item **Magnet & Cooling - Magnet Status**.

Fig. 160: Checking the Ceramic Carbon Resistors



Tab. 9 Expected CCR resistor values

| 41-pin plug, measured in Ohms       | Expected values<br>room temperature | Expected values<br>77 K temperature | Expected values<br>4.2 K temperature |
|-------------------------------------|-------------------------------------|-------------------------------------|--------------------------------------|
| All ceramic carbon resistors (CCRs) | 880 - 1050                          | 1150 - 1350                         | 2800 - 3380                          |





If the host (or laptop), is not available, check CCR values directly from the pins (a-e) on the 41-pin connector on the turret. (→ TSG - Checking the Carbon Resistors / Page 160)



The magnet helium level can be recovered if the temperature of the CCR is still colder than nitrogen temperature.

The CCR (CARBON RES 1) reading must be in the range of between 1200 and 3400 Ohms.



If temperatures are between N2 and the helium temperature, a BOTTOM FILL must be performed IMMEDIATELY.

Fig. 161: OR98 Magnet Acceptance Document (example)

| <b>SIEMENS</b>                                                                                                                          |                         | <b>Magnet Acceptance Documentation</b> |                      |
|-----------------------------------------------------------------------------------------------------------------------------------------|-------------------------|----------------------------------------|----------------------|
|                                                                                                                                         |                         | <b>Installation Data</b>               |                      |
| <b>MR Magnet Technology Quality Regulations</b>                                                                                         |                         |                                        |                      |
| <b>Cryostat Number</b>                                                                                                                  | 21236477                |                                        | <b>SMT Reference</b> |
| 21245226 05.10                                                                                                                          |                         |                                        |                      |
| <b>System Description</b>                                                                                                               |                         |                                        |                      |
| 1.495T Actively Shielded Superconducting Twin Turret Magnet System                                                                      |                         |                                        | OR98                 |
| <b>Measured Magnet Centre (Offset)</b>                                                                                                  |                         |                                        | +1 (P) mm            |
| <b>Gradient Coil Position</b>                                                                                                           |                         |                                        | 39 mm                |
| <b>Overall System Length</b>                                                                                                            |                         |                                        | 1409 mm              |
| This system requires 12 Z2 rings and 2.38 Kg of iron to shim. This is a theoretical calculation that may change due to site conditions. |                         |                                        |                      |
| <b>Operating Parameters</b>                                                                                                             |                         |                                        |                      |
| <b>Description</b>                                                                                                                      |                         | <b>Measured value</b>                  |                      |
| Operating Current @ 1.495T                                                                                                              |                         | 510.83 AMPS                            |                      |
| Field/Current Ratio                                                                                                                     |                         | 12.46 kHz/0.1 Amp                      |                      |
| Magnet Type: OR98-??                                                                                                                    |                         |                                        |                      |
| <b>Carbon Resistor Values</b>                                                                                                           |                         |                                        |                      |
| <b>POSITION</b>                                                                                                                         | <b>Room Temperature</b> | <b>Helium Temperature</b>              |                      |
| CCR4 TDC Outer (D-E)                                                                                                                    | 926                     | 2951                                   |                      |
| CCR3 Top Inner (C-E)                                                                                                                    | 1060                    | 2812                                   |                      |
| CCR2 Middle Inner (B-E)                                                                                                                 | 1035                    | 2977                                   |                      |
| CCR1 Bottom Inner (A-E)                                                                                                                 | 1093                    | 3138                                   |                      |
| TDC Calibration Values:                                                                                                                 | 882.4                   | 1181.8                                 |                      |
| <b>Helium Level Probe Resistance (When using a 250 mA supply)</b>                                                                       |                         |                                        |                      |
|                                                                                                                                         | Probe 1                 | Probe 2                                |                      |
| 0% Resistance                                                                                                                           | 107.0                   | 107.0                                  | Ohms                 |
| 0 % Calibration                                                                                                                         | 448                     | 448                                    |                      |
| 100 % Resistance                                                                                                                        | 17.9                    | 17.9                                   | Ohms                 |
| 100 % Calibration                                                                                                                       | 75                      | 75                                     |                      |
| Diode Table No                                                                                                                          | -2                      |                                        |                      |
| Pressure Heater Res                                                                                                                     | 22.1                    | Ohms                                   |                      |
| <b>Final Temperature Sensor Recon Values</b>                                                                                            |                         |                                        |                      |
| 50Kbore                                                                                                                                 | 50K Link                | Coldhead                               | Recon Margin         |
| 56                                                                                                                                      | 55                      | 41                                     | 782 Milliwatts       |
| Checked by                                                                                                                              | K Suraina               |                                        | Signature            |
| Date                                                                                                                                    | 28/05/2010              |                                        |                      |

Fig. 162: OR99 Magnet Acceptance Document (example)

|                                                                                                                                         |                            |                                                              |                      |                |  |
|-----------------------------------------------------------------------------------------------------------------------------------------|----------------------------|--------------------------------------------------------------|----------------------|----------------|--|
| <b>SIEMENS</b>                                                                                                                          |                            | <b>Magnet Acceptance Documentation<br/>Installation Data</b> |                      |                |  |
| <b>MR Magnet Technology Quality Regulations</b>                                                                                         |                            |                                                              |                      |                |  |
| <b>Cryostat Number</b>                                                                                                                  | 21240447                   |                                                              | <b>SMT Reference</b> | 21243544 05.10 |  |
| <b>System Description</b>                                                                                                               |                            |                                                              |                      |                |  |
| 2.894T Actively Shielded Superconducting Twin Turret Magnet System                                                                      |                            |                                                              |                      | OR99           |  |
| Measured Magnet Centre (Offset)                                                                                                         |                            |                                                              | -3 (S)               | mm             |  |
| Gradient Coil Position                                                                                                                  |                            |                                                              | 43                   | mm             |  |
| Overall System Length                                                                                                                   |                            |                                                              | 1669 mm              |                |  |
| This system requires 12 Z2 rings and 1.44 Kg of iron to shim. This is a theoretical calculation that may change due to site conditions. |                            |                                                              |                      |                |  |
| <b>Operating Parameters</b>                                                                                                             |                            |                                                              |                      |                |  |
| <u>Description</u>                                                                                                                      | <u>Specification value</u> | <u>Measured value</u>                                        |                      |                |  |
| Operating Current @ 2.894T                                                                                                              | < 500 AMPS                 | 465.97 AMPS                                                  |                      |                |  |
| Field/Current Ratio                                                                                                                     | N/A                        | 26.44 kHz/0.1 Amp                                            |                      |                |  |
| Magnet Type: OR99-??                                                                                                                    |                            |                                                              |                      |                |  |
| <b>Carbon Resistor Values</b>                                                                                                           |                            |                                                              |                      |                |  |
| <u>POSITION</u>                                                                                                                         | <u>Room Temperature</u>    | <u>Helium Temperature</u>                                    |                      |                |  |
| CCR4 TDC (D-E)                                                                                                                          | 1072                       | 2863                                                         |                      |                |  |
| CCR3 Top (C-E)                                                                                                                          | 872                        | 3114                                                         |                      |                |  |
| CCR2 Middle (B-E)                                                                                                                       | 1027                       | 3182                                                         |                      |                |  |
| CCR1 Bottom (A-E)                                                                                                                       | 1036                       | 2963                                                         |                      |                |  |
| TDC Calibration Values:                                                                                                                 | 1022.6                     | 1322.1                                                       | 2991.4               |                |  |
| <b>Helium Level Probe Resistance (When using a 250 mA supply)</b>                                                                       |                            |                                                              |                      |                |  |
|                                                                                                                                         | Probe 1                    | Probe 2                                                      |                      |                |  |
| 0% Resistance                                                                                                                           | 80.3                       | 80.3                                                         | Ohms                 |                |  |
| 0 % Calibration                                                                                                                         | 336                        | 336                                                          |                      |                |  |
| 100 % Resistance                                                                                                                        | 0.5                        | 0.5                                                          | Ohms                 |                |  |
| 100 % Calibration                                                                                                                       | 002                        | 002                                                          |                      |                |  |
| Diode Table No                                                                                                                          | -2                         |                                                              |                      |                |  |
| Pressure Heater Res                                                                                                                     | 20.5                       | Ohms                                                         |                      |                |  |
| <b>Final Temperature Sensor Recon Values</b>                                                                                            |                            |                                                              |                      |                |  |
| 50Kbore                                                                                                                                 | 50K Link                   | Coldhead                                                     | Recon Margin         |                |  |
| 54                                                                                                                                      | 50                         | 40                                                           | K 874 Milliwatts     |                |  |
| Checked by                                                                                                                              | K Suraina                  |                                                              | Signature            |                |  |
| Date                                                                                                                                    | 10/05/2010                 |                                                              |                      |                |  |
| Page 1 of 1                                                                                                                             |                            |                                                              | 25-AD EPS 000 01     |                |  |

Fig. 163: OR103 Magnet Acceptance Document

| SIEMENS                                                                    |                                                                          | Magnet Acceptance Documentation          |                           |
|----------------------------------------------------------------------------|--------------------------------------------------------------------------|------------------------------------------|---------------------------|
|                                                                            |                                                                          | Installation data                        |                           |
|                                                                            |                                                                          | MR Magnet Technology Quality Regulations |                           |
| <b>Cryostat number</b>                                                     | 33002284                                                                 | <b>SMT Shipped System</b>                | 33002281 07.11            |
| <b>System description</b>                                                  | 2.894T Actively Shielded Superconducting Twin Turret Magnet System OR103 |                                          |                           |
|                                                                            | <b>SKYRA</b>                                                             | <b>MONET</b>                             |                           |
| <b>Measured magnet centre (Offset)</b>                                     | -4 (S)                                                                   | -4 (S)                                   | mm                        |
| <b>Gradient Coil Position</b>                                              | 174                                                                      | 174                                      | mm                        |
| <b>Overall Magnet Length</b>                                               | 1670                                                                     | 1670                                     | mm                        |
| <b>This system requires</b>                                                | 0                                                                        | 0                                        | Z2 Rings                  |
| <b>This system requires</b>                                                | 2.63                                                                     | 2.63                                     | kg of iron to shim        |
| (This is a theoretical calculation that may change due to site conditions) |                                                                          |                                          |                           |
| <b>Operating parameters</b>                                                |                                                                          |                                          |                           |
| <b>Description</b>                                                         | <b>Specification value</b>                                               | <b>Measured value</b>                    |                           |
| Operating current @ 2.894T                                                 | < 500 Amps                                                               | 487.35 Amps                              |                           |
| Field/Current Ratio                                                        | N/A                                                                      | 25.28 kHz/0.1 Amp                        |                           |
| Magnet type:                                                               | OR103-00                                                                 |                                          |                           |
| <b>Carbon resistor values</b>                                              |                                                                          |                                          |                           |
| <b>Position</b>                                                            | <b>Room temperature</b>                                                  | <b>Nitrogen temperature</b>              | <b>Helium temperature</b> |
| CCR4 d - e top                                                             | 1071 Ohms                                                                |                                          | 2573 Ohms                 |
| CCR3 c - e                                                                 | 991 Ohms                                                                 |                                          | 2612 Ohms                 |
| CCR2 b - e                                                                 | 929 Ohms                                                                 |                                          | 2797 Ohms                 |
| CCR1 a - e bottom                                                          | 1007 Ohms                                                                |                                          | 2841 Ohms                 |
| CCR4 calibration values                                                    | 1031 Ohms                                                                | 1304 Ohms                                | 2621 Ohms                 |
| <b>Helium level probe resistance (when using a 250mA supply)</b>           |                                                                          |                                          |                           |
|                                                                            | <b>Probe 1</b>                                                           | <b>Probe 2</b>                           |                           |
| 0% resistance                                                              | 107 Ohms                                                                 | 107 Ohms                                 |                           |
| 0% calibration                                                             | 448                                                                      | 448                                      |                           |
| 100% resistance                                                            | 10.8 Ohms                                                                | 10.8 Ohms                                |                           |
| 100% calibration                                                           | 45                                                                       | 45                                       |                           |
| Diode table no.                                                            | 2                                                                        |                                          |                           |
| Pressure heater resistance                                                 | 20.5 Ohms                                                                |                                          |                           |
| <b>Final temperature sensor recon values</b>                               |                                                                          |                                          |                           |
| Shield bore                                                                | Shield Link                                                              | Coldhead                                 | Recon margin              |
| 55K                                                                        | 48K                                                                      | 39K                                      | 853mW                     |
| <b>Checked by</b>                                                          | <b>Signed</b>                                                            |                                          |                           |
| K. Baker                                                                   | <i>Km Baker</i>                                                          |                                          |                           |
| <b>Date</b> 11/07/2011                                                     |                                                                          |                                          |                           |
| Page 1 of 1                                                                |                                                                          | 25-AD EPS 000 01                         |                           |




#### 4.2.3.6 Check the Magnet Pressure

Check that the helium vessel has a positive pressure:

- Check for positive pressure on the pressure gauge.
- Check the pressure at the MRC (MRAWP) Host.


Fig. 164: MSUP - Pressure heater LED

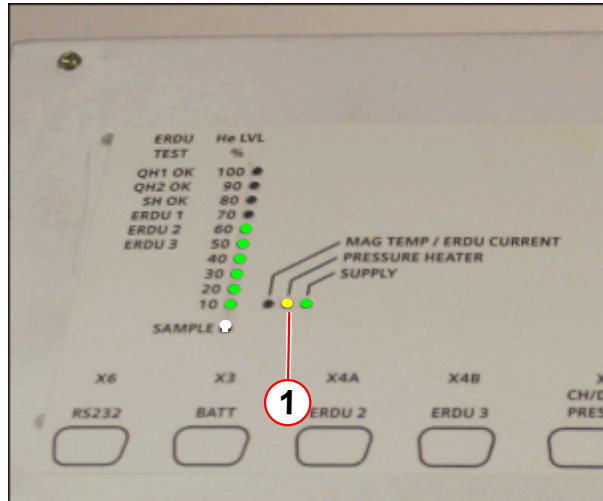

(1) Pressure heater LED - yellow


- Check to ensure that the PRESSURE HEATER LED on the MSUP lights up every second to maintain a positive pressure.
  - If the LED is always on then it is likely that a burst has blown or there is a gross leak.
  - If the LED does not come on then the pressure is too high or the pressure heater circuit is faulty.
- Access the Service Software at the MRC (MRAWP) Host for more information.
- If the pressure in the magnet falls too low, the coldhead will switch off as it approaches atmospheric pressure.

Fig. 165: OR98 - Pressure gauge

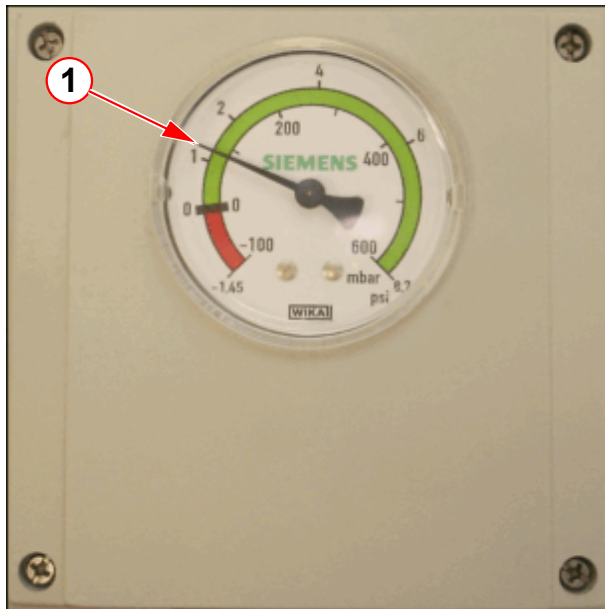

(1) Expected gauge pressure 0.7 - 2.3 psi (G)

Fig. 166: OR103 - Pressure gauge (expected pressure)

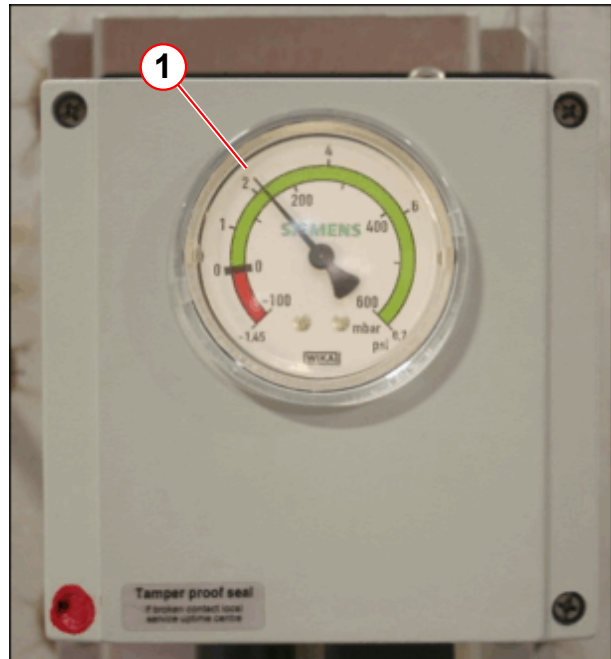


(1) Expected gauge pressure 1.7 - 3.3 psi (G)

- If the pressure gauge is within the **RED zone**, the magnet is sub-atmospheric.
  - » **DO NOT** open the magnet bypass valve, or disturb any high pressure fitting.

#### 4.2.3.7 Deactivate the Pressure Heater

- Access the Service Software at the MRC (MRAWP) Host.
- Select the menu item Magnet & Cooling - Magnet Status.
- At MSUP Control - Set the pressure heater mode to **OFF**.
- Select **SAVE**.

Fig. 167: Service screen Magnet & Cooling - Magnet status - MSUP Control

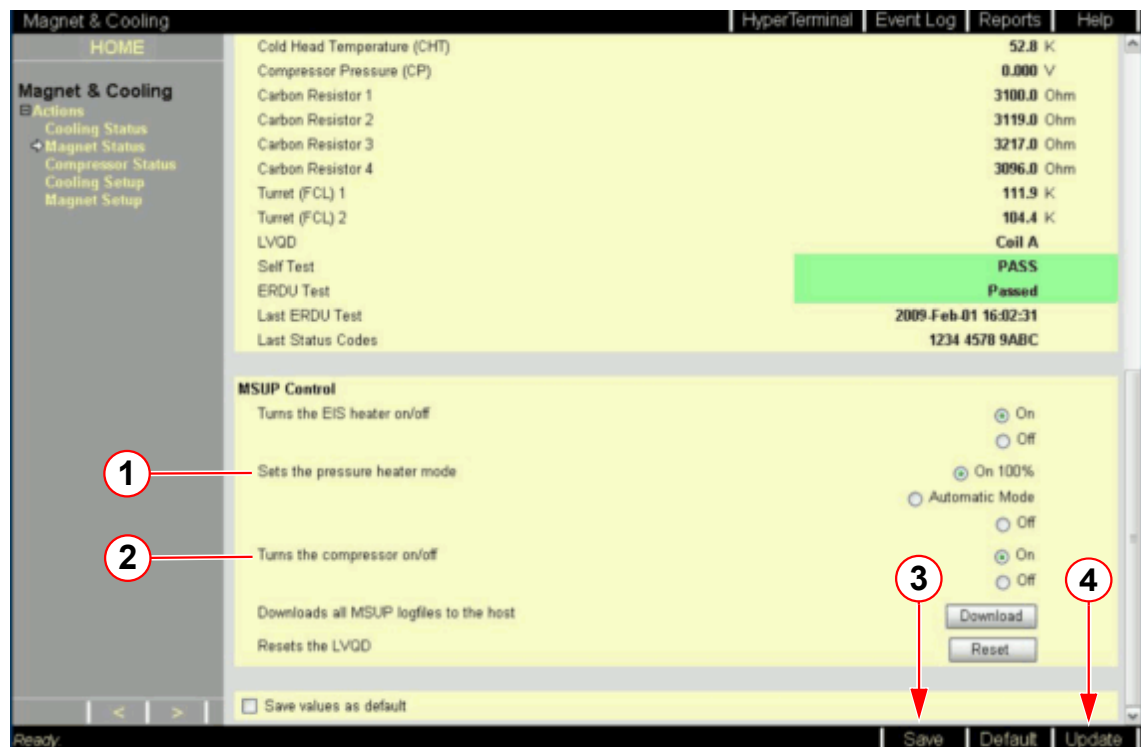

- (1) Sets the pressure heater mode
- (2) Sets the compressor on/off
- (3) Save button
- (4) Update button



The pressure heater will reactivate within 60 min.



## 4.2.3.8 Depressurize the Helium Vessel


**CAUTION**

**Risk of ice formation.**

The magnet bypass valve is only **OPENED** for depressurization of the helium vessel. The bypass valve is **CLOSED** during normal operation.

- » To avoid air ingress, **CLOSE** the magnet bypass valve immediately after the helium vessel is depressurized.


Fig. 168: Bypass valve positions for depressurization

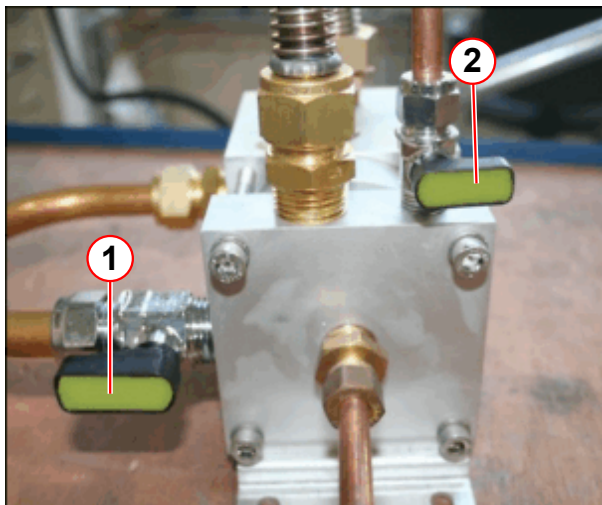

- (1) Magnet bypass valve OPEN
- (2) Cold head bypass valve CLOSED

- Fully open the magnet bypass valve to depressurize the helium vessel.
- Wait until zero psi is displayed on the pressure gauge.
- The pressure switch will activate.
  - » the cold head will turn off.
- **Close** the magnet bypass valve to prevent air ingress.
- Warm up the vent valve assembly.
  - » Take care not to damage any components when using a heat gun.

Fig. 169: MSUP - Magnet pressure is low

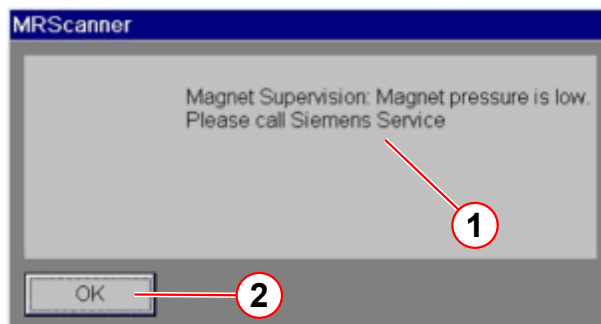

- (1) Pop-up message
- (2) Click OK to clear

On the MR Scanner:

- A pop-up alarm message is displayed.
  - » Magnet pressure is low.
  - » Click OK to clear the alarm.
- The alarm box will sound when pressure falls to alarm limits.
  - » The yellow WARNING LED lights on the alarm box.
  - » The alarm buzzer is activated
  - » Press AUDIO ALARM OFF to silence the alarm.


## 4.2.4 Checking the Quench Valve and Burst Disk

### 4.2.4.1 Quench Valve Leak Check

Post quench if there is no magnet pressure a burst disk has ruptured. Check the graphite (quench valve) burst disk.

- Warm up the quench valve and venting to room temperature, using a heat gun.
- The elbow assembly must be removed for inspection and leak check of the burst disk.

Fig. 170: Horizontal service port

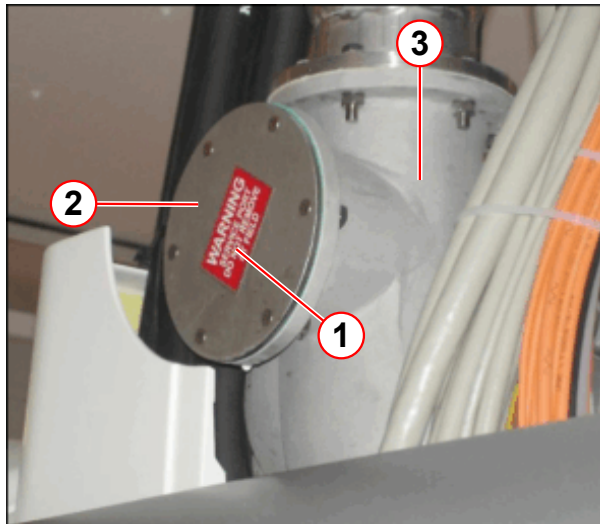

- (1) Warning label
- (2) Port (blanked)
- (3) Quench elbow


#### OR98 only - Horizontal Service Port

- This port cover must NOT be removed.
- **Do not** use it as an inspection cover to check the quench valve burst disk.

It may be necessary to remove the complete elbow from the acoustic de-coupler. This will depend on the system installation and how cables are routed. It may also be necessary to remove:

- The vent line from the magnet bypass valve.
- The vent line from the OVC burst disk port (if fitted).

Fig. 171: Quench valve removal

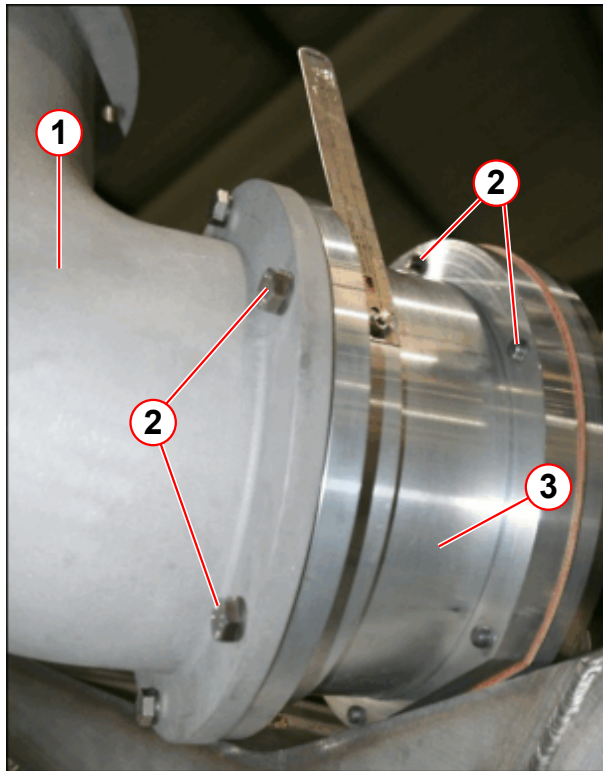


- (1) Elbow
- (2) M8 screws and nordlock washers (6 off) each end
- (3) Quench valve assembly

The elbow is secured to the quench valve assembly with 6 x M8 screws and nordlock washers.

- Disconnect the quench valve from the elbow:
  - Remove the (6 off) M8 screws and washers from the elbow.


Fig. 172: Check the quench valve burst disk

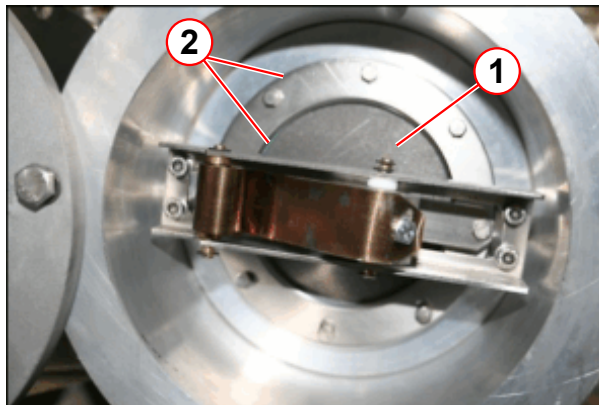

- (1) Graphite burst disk
- (2) Leak check retaining ring

- Inspect the graphite burst disk for breakage, or cracks.
- If the burst disk has ruptured, or is cracked, it must be replaced **immediately**.
- Carefully warm and dry around the burst disk retaining ring and spring mechanism ready for leak checking.

Fig. 173: Water corrosion (example)

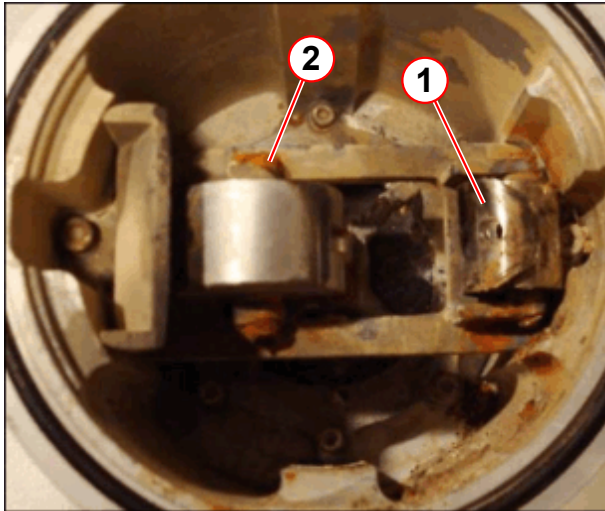

- (1) Spring broken  
(2) Corroded screws

Check for ice and water corrosion or breakage of:

- The spring plate assembly.
- All screws.
- Corroded quench valve components must be replaced.
- Check that the quench valve drain hole is clear of debris and that the drain screw is removed.

Complete leak check:

- With no damage to the burst disk or spring plate visible, continue to leak check.
- Check that the magnet pressure is a minimum of 0.5 psi (G) for leak checking.
- Starting at the bottom of the disk, slowly leak check around the retaining ring. If there is no leak found the O-ring behind the spring plate is OK.


#### 4.2.4.2 Quench Valve O-ring Inspection

- To inspect a leaking or damaged O-ring the burst disk spring plate must be removed from the quench valve.
  - » Check that there is no pressure in the magnet.
  - » Depressurize if necessary.



The spring exerts significant force on these screws.


Fig. 174: Quench valve burst disk assembly

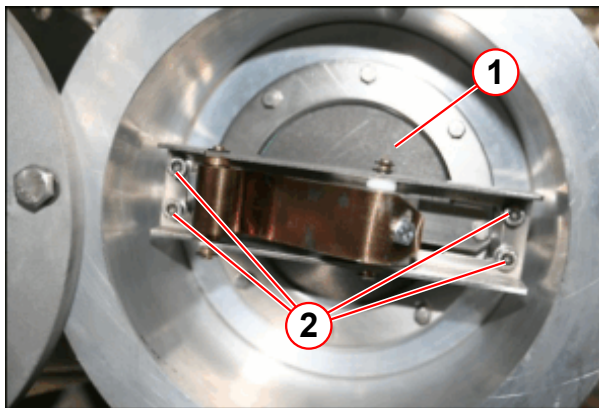

- (1) Graphite burst disk
- (2) Spring plate retaining screws M5 x 4

- Diametrically unscrew and remove the (4 off) M5 cap head screws from the spring retaining plate.
- Carefully lift out the burst disk spring plate assembly from the quench valve body.


Fig. 175: Quench valve burst disk O-ring

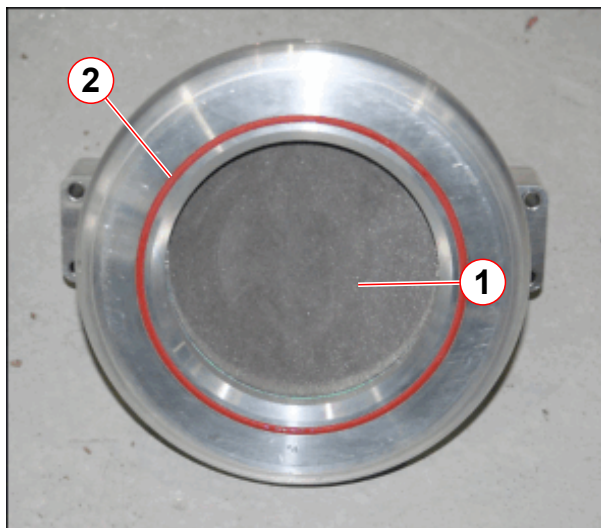

- (1) Burst disk
- (2) O-ring

- Inspect the quench valve O-ring.
  - » Replace if damaged.
- Clean the O-ring groove and O-ring before re-fitting.
- Check that the O-ring is seated evenly all around the groove.
  - » **DO NOT** use any form of grease on the O-ring.

Fig. 176: Check the quench valve burst disk

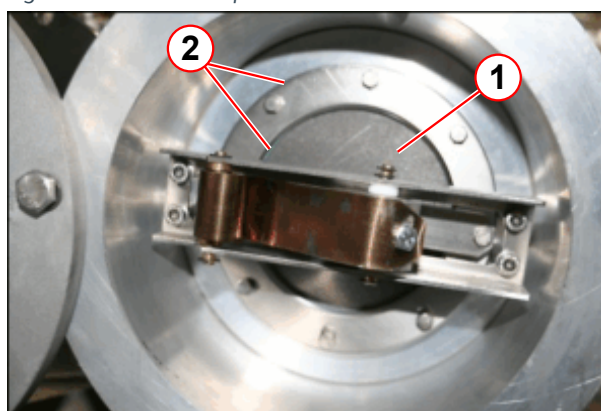


- (1) Graphite burst disk
- (2) Leak check retaining ring

- Re-fit the burst disk assembly into the quench valve body with the 4 x M5 screws.
- Check that the burst disk remains level with the sealing face.
- Torque 4 x M5 screws in diametric pattern 6 Nm.
- Check that the magnet pressure is a minimum of 0.5 psi (G) for leak checking.
- Leak check around the retaining ring.



Fig. 177: Quench valve and elbow

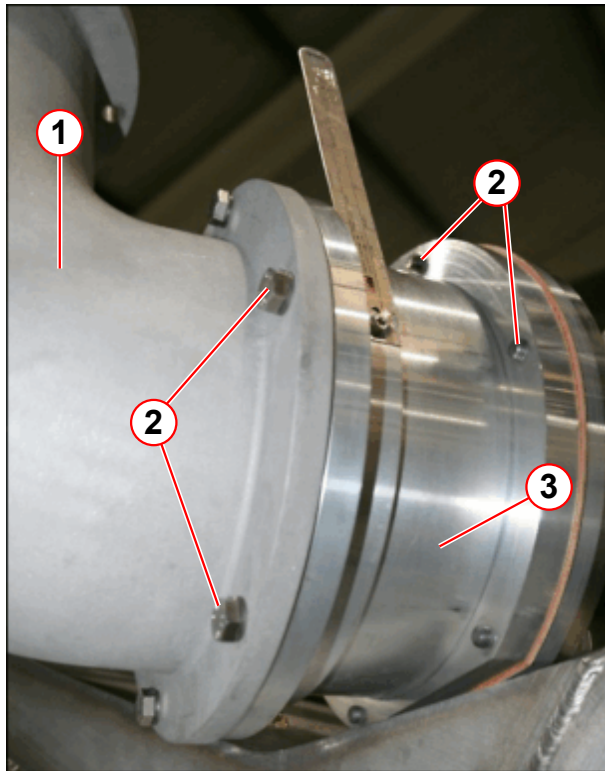

- (1) Elbow
- (2) M8 screws and nordlock washers (6 off) each end
- (3) Quench valve assembly

If all the leak check and burst disk inspections are complete:


- Refit the elbow is to the quench valve assembly with 6 x M8 screws and nordlock washers.
- If removed from the acoustic decoupler, refit the acoustic decoupler.

#### 4.2.5 Checking the Auxiliary Burst Disk

- Check that the magnet is at zero pressure.
- If the elbow has not been removed for the quench valve inspection:
  - Warm up the auxiliary burst disk assembly to room temperature, using a heat gun.



Fig. 178: Auxiliary burst disk assembly

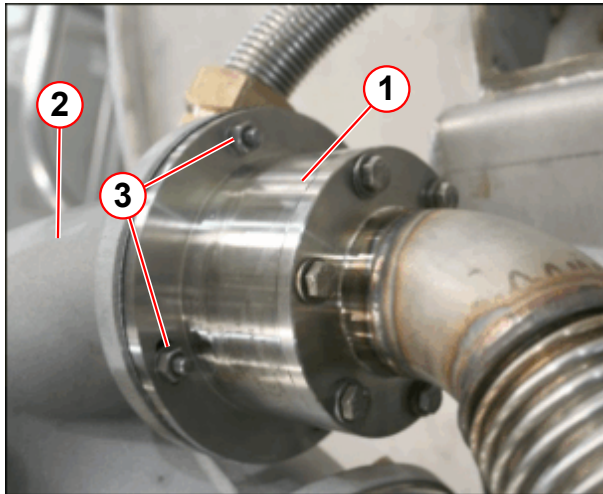


- (1) Auxiliary burst disk
- (2) Elbow
- (3) M6 screws, washers and nuts (4 off)

- Disconnect the auxiliary burst disk from the elbow
  - Remove the M6 screws, washers and nuts (4 off).
- Remove the vent pipe from the bypass assembly at the quench elbow.
- **Immediately** place a ball of paper towelling into the open elbow port.



Do not touch the metal membrane; it is very fragile.


Fig. 179: Auxiliary burst disk

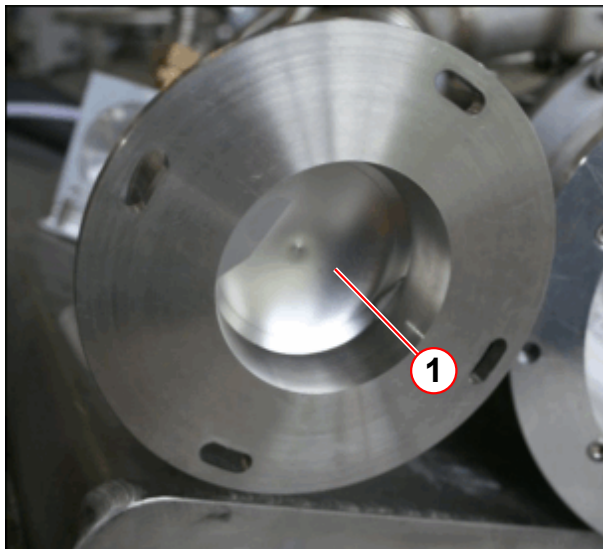

- (1) Burst disk

- Visually inspect the auxiliary burst disk membrane to ensure that it is not damaged.
- A damaged or ruptured burst disk must be replaced immediately. If a new burst disk is not available the vent line must be blanked until a replacement can be obtained.
  - » Fit the 2 x red blanking plates or the gasket blanking seal (from spares kit).

Fig. 180: Auxiliary burst disk gasket

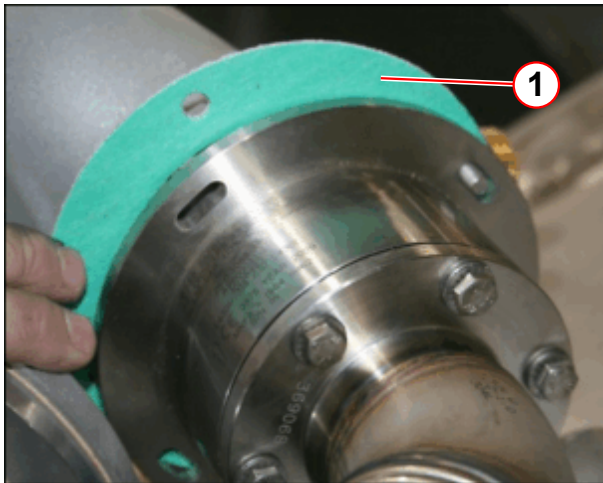

(1) Gasket

If no damage is seen to the burst disk:

- Ensure that the gasket seal is not damaged, and that the surface is clean.
- Fit the auxiliary burst disk assembly to the elbow.
  - use the M6 screws, washers and nuts.
  - torque the M6 screws to 6 Nm.

#### 4.2.5.1 Blanking the Auxiliary Burst Disk Vent Line

In the event that the auxiliary burst disk has ruptured during a quench the vent path must be sealed from air ingress. If there is no replacement burst disk available the auxiliary vent must be sealed with a blanking gasket.

There are 2 options available for blanking the auxiliary burst disk and vent line. (Dependent on system status).

1. 2 x red blanking plates.
2. 1 x gasket blanking seal.



#### WARNING

##### Safety device removed

The magnet must not be ramped until a new auxiliary burst disk is fitted.

- » Do not ramp with the auxiliary blanking gasket fitted.
- » Do not ramp with a ruptured auxiliary burst disk.



#### WARNING


##### Risk of air ingress

If the auxiliary burst disk is broken, IMMEDIATELY fit the blanking gasket seal.

- » The blanking gasket seal is supplied in the magnet spares kit.
- » Order a replacement auxiliary burst disk and fit before ramping the magnet.

Red blanking plate fitment

Fig. 181: Red Blanking plate fitted

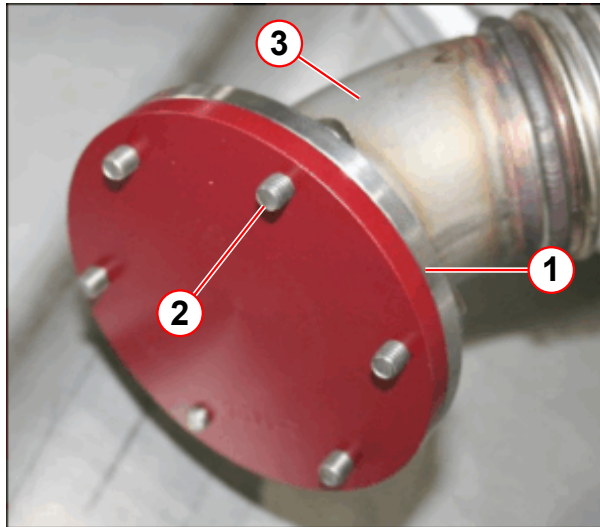

- (1) Leak check seal
- (2) M5 screws and washers (6 off)
- (3) Auxiliary vent line

- 2 x red anodized blanking plates and O-rings are provided in the magnet site box.
- If the burst disk is damaged, remove and fit the red blanking plates immediately to
  - The auxiliary vent line.
  - The elbow.
- A replacement burst disk and seal must be obtained urgently.
- With the blank fitted allow magnet to pressurize.
- Leak check red blanking plate on the end of the auxiliary vent line (Pos.1).


### Gasket seal blank

Fig. 182: Gasket blank and label



- Check that the warning label is secured to the gasket seal with a cable tie.

Fig. 183: Blanking gasket fitted

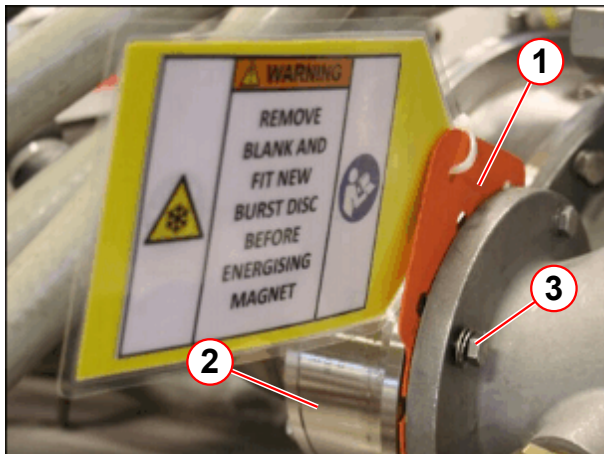

- (1) Blanking gasket seal
- (2) Auxiliary burst disk (ruptured)
- (3) M6 screws x 4

Leave the ruptured burst disk attached to the auxiliary turret vent line.

- Remove the 4 x M6 screws, washers and nuts retaining the burst disk to the elbow.
- Remove the gasket seal.
- Fit the blanking gasket seal between the burst disk and the elbow.
- Fit the auxiliary burst disk and gasket blanking seal to the elbow.
  - Use the M6 screws, washers and nuts.
  - Torque the M6 screws to 6 Nm.

When the blanking gasket is fitted:

- Warm up the auxiliary venting using a heat gun
- Allow the magnet to pressurize.
- Leak check around the joint on the high pressure turret side.



When a replacement auxiliary burst disk is available, remove the blanking gasket. Fit the new burst disk following the **Replacement procedures**.

## 4.2.6 Final Work Steps

### 4.2.6.1 Reactivate the Pressure Heater

- Access the Service Software at the MRC (MRAWP) Host.
- Select the menu item Magnet & Cooling - Magnet Status.
- At MSUP Control - Set the pressure heater mode to **On 100%**. (→ Fig. 47 Page 46)
- Select SAVE.
- Monitor the magnet pressure.
  - Select Update to refresh data until pressure reaches  $\geq 15.5$  psi (A).



Data is only accurate from the last refresh/update of the page, or from when the page was opened last.




During pressurization, the cold head will switch on automatically. If the cold head was turned off at MRC, select 'pressure heater status'.

- At MSUP Control - Set the pressure heater mode to Automatic mode.

- Select SAVE.

Fig. 184: Check magnet pressure level

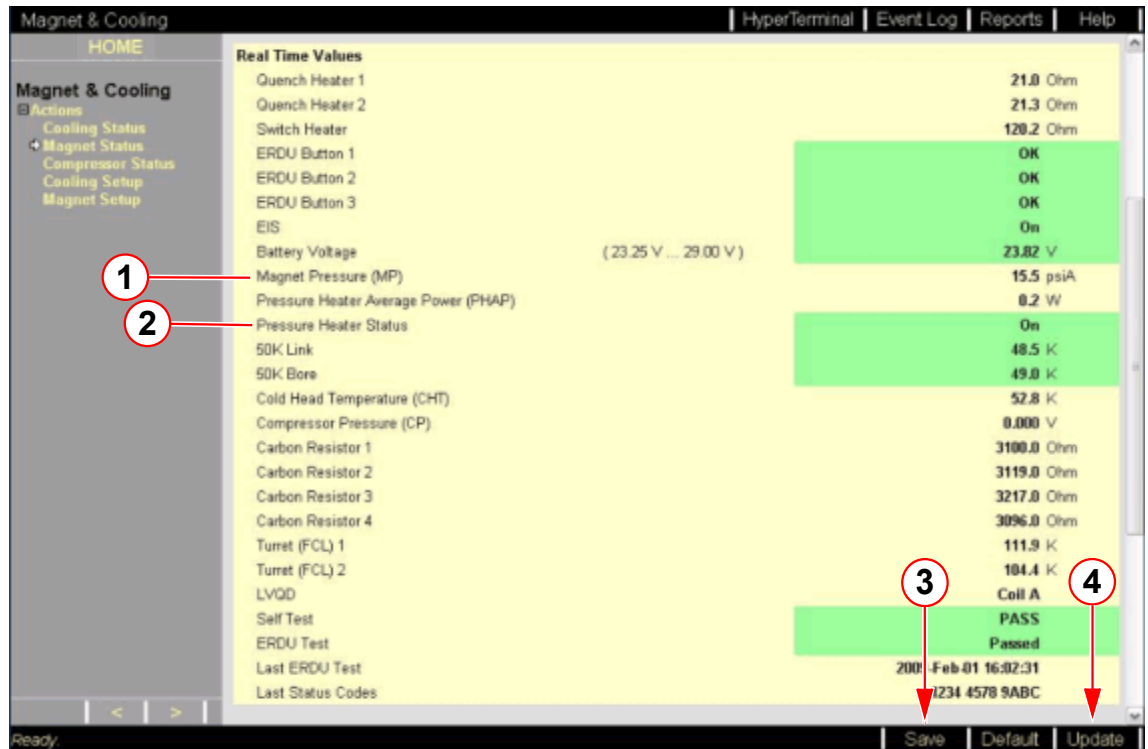


- (1) Magnet pressure (MP)
- (2) Pressure heater status
- (3) Save button
- (4) Update button

#### 4.2.6.2 Final Checks

- Helium fill to the required helium level for ramping.
- Order replacement components if the spare parts have been used (burst disk, seals).
- Perform a leak check of the turret venting and components.
  - » All vent pipes and flanges on the high pressure side of the helium vessel.



Fig. 185: Turret leak check

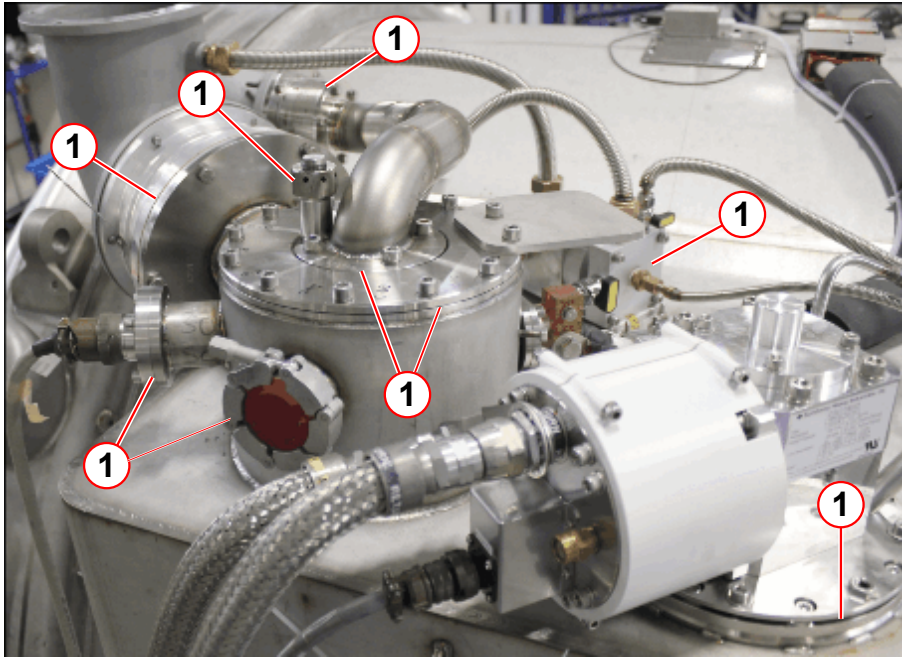

(1) Leak check points



Allow a minimum of 2 hours after filling before ramping the magnet.

#### 4.2.6.3 Quench Report

- Fill out the Quench report. - MR - Quality Database and Quench Report.

#### 4.2.6.4 TSG - Checking the Carbon Resistors

If required, the Ceramic Carbon Resistors (CCR) values can be checked from the MSUP or manually by using a multimeter.

- The measurements are taken across pins **a-e**, **b-e**, **c-e**, **d-e** on the 41-pin turret connector.



### CAUTION

**Risk of asphyxiation.**

**Failure to adhere to safety precautions may result in serious bodily injury.**

» **The magnet must be RAMPED DOWN before attempting this procedure.**

#### Tooling requirements

- Multimeter and pin leads.
- Laptop/PC with Terminal Emulation software, (PuTTY).




- RS232 cable.

### Procedure

Using a multimeter:


- Connect the pin leads.
- Turn the mode selector to resistor mode  $\Omega$ .

*Fig. 186: Turret - 41pin connector*

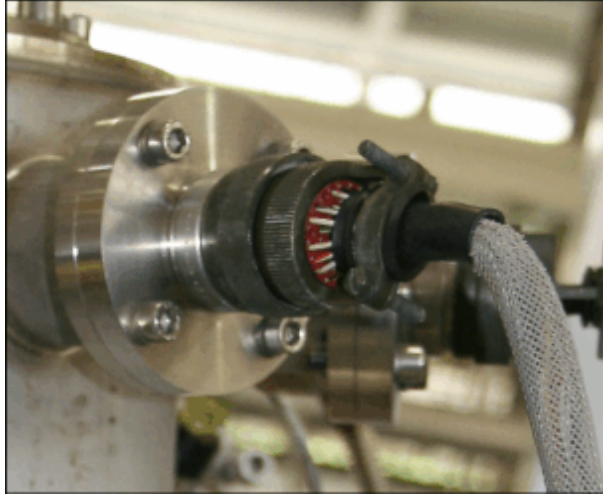

- Disconnect the 41-pin connector.
- Measure across pins **a-e**, **b-e**, **c-e**, **d-e**
- Compare the values with the Magnet Acceptance Document.
- Re-connect the 41-pin connector.

Fig. 187: 41pin turret connector - ceramic carbon resistor pins (room temperature values)

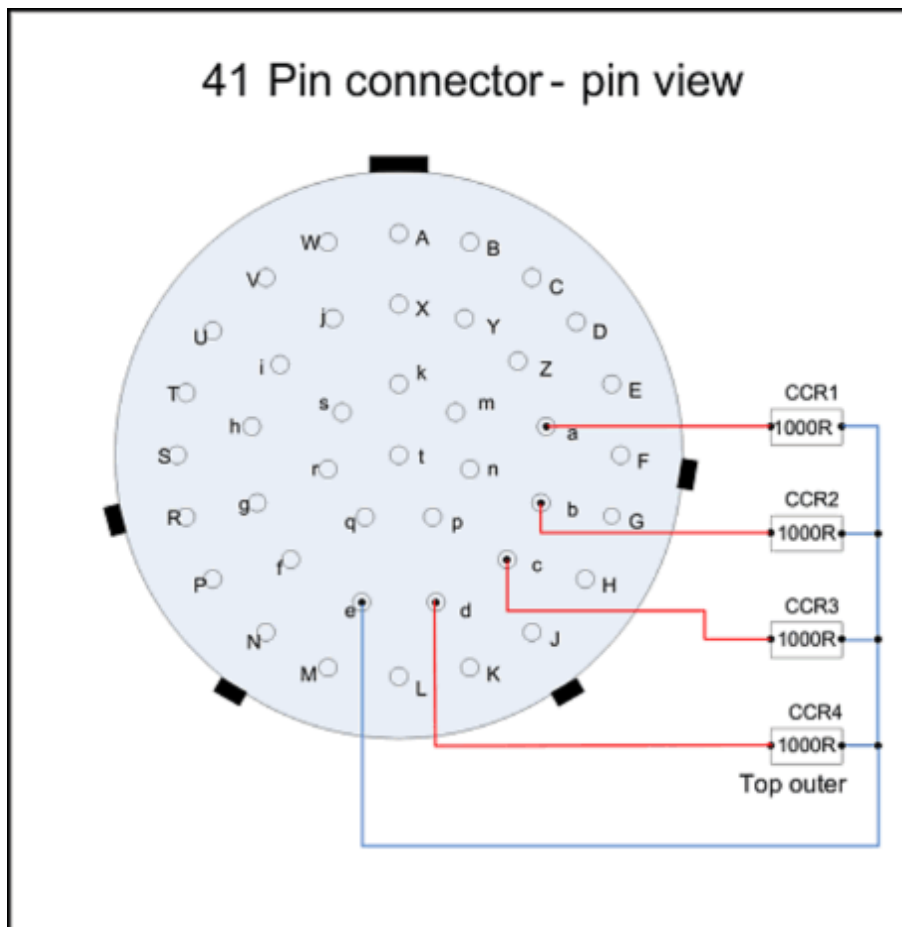

## 5.1 Prerequisites for Shimming

### 5.1.1 Only for Systems with a Tim Dockable Patient Table



If the system has only one dockable patient table, no further action is necessary.

If the system has more than one dockable patient table, you have to check which patient table is registered.

Only one dockable patient table is registered in the factory.

A correctly registered patient table is necessary for shimming.

#### 1. Check which dockable patient table is already registered.

- Open a command tool and enter the command:

```
C:\Users\meduser>MrVerifyPTAB
```

#### Example:

Fig. 188: MrVerifyPTAB

```
86); Numaris 4 VD13A ; 25-Sep-2011 00:51; .
```

```
-----
current system settings:
```

```
mobile.....:1
frequency.....:123200000Hz
field strength.....:2.9T|
```

```
current PTAB settings read from HW:
```

```
material number.....:10433522
serial number.....:1002
```

```
PTAB devices registered on host:
```

```
device 0:
material number.....:10433522
serial number.....:1002
-----
```

```
finished successfully.
```

#### 2. Use the registered dockable patient table for shimming.



Check the correct adjustment of the dockable patient table.

Refer to the Installation Instructions (→ Tim Dockable Table (Mobile Patient Table) / M7-010.812.01) if an adjustment is necessary.

If the system has a 2nd or more dockable patient tables, these additional patient tables need to be registered in the system software.

For detailed information, please refer to document (→ Additional Tim Dockable Table / MR07-000.814.05)


## 5.2 Passive Shim



Information regarding passive Shim is not published as hardcopy documentation. For more information, refer to the Online Help

- Start the Shim workflow from the SeSo platform.

Fig. 189: Service Software

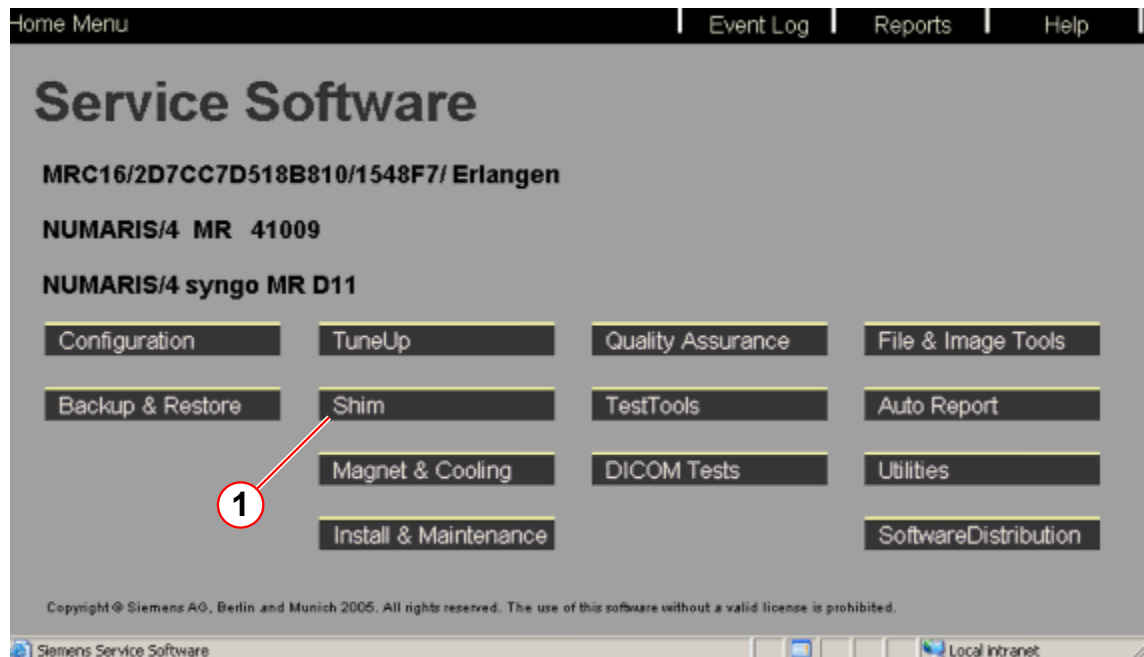

(1) Shim workflow

## 6.1 Changes to previous version

### Version 02

- Chapter 03
  - Renaming of the chapter.
  - Update of the chapter.
  - Note of the Quench line acceptance protocol added.
- Chapter 04
  - Note of the Quench line acceptance protocol added.

### Version 03

- Chapter 03
  - Update of links
- Chapter 04
  - Update of links
  - Update of MRC Host connection.

### Version 04

- Chapter 03
  - Venting Leak Check. Reference to use 5 Nm added.
- Chapter 04
  - What to do after Quench. Mobile reference deleted.

### Version 05

- Chapter 03
  - New leak detector added.
- Chapter 05
  - Dockable patient table added.

### Version 06

- Chapter 04
  - MNO Frequency range added.
  - Update of the screenshots

### Version 07

- Chapter 04
  - Update of note: functionality of host computer.
  - Update of terminal emulation software.

### Version 08

- All,
  - Editorial changes
- Chapter 03, Perquisites Ramping:
  - New path to the acceptance protocol and deleted the sign off by the project manager, Charm MR\_00455599)



- Chapter 04, Ramping the Magnet:
  - Transport test, updated. Charm MR\_00455599)
  - Warning for visual inspection of quench tube added. (→ OR98/99/103 Ramping / Page 93)

**Version 09**

- Chapter 04, Ramping the Magnet:
  - Co-Siting of Magnets added.

**Version 10**

- Chapter 03, Prerequisites for Ramping the Magnet
  - Update of "Venting leak check"
- Chapter 04, Ramping the Magnet:
  - Update of "What to do after a quench" (ID 201048, ID 212345)

There are no Hazard IDs in this document.

All documents may only be used for rendering services on Siemens Healthcare Products. Any document in electronic form may be printed once. Copy and distribution of electronic documents and hardcopies is prohibited. Offenders will be liable for damages. All other rights are reserved.

[healthcare.siemens.com/services](http://healthcare.siemens.com/services)

Siemens Healthineers Headquarters  
Siemens Healthcare GmbH  
Henkestr. 127  
91052 Erlangen  
Germany  
Telephone: +49 9131 84-0  
[siemens.com/healthineers](http://siemens.com/healthineers)

Print No.: M7-010.815.20.10.02 | Replaces: M7-010.815.20.09.02  
Doc. Gen. Date: 10.17 | Language: English  
© Siemens Healthcare GmbH, 2010

**[siemens.com/healthineers](http://siemens.com/healthineers)**