

820-0241 – Filter Replacement for Cryocooler Compressor Valves P/N: 205-0121-13 (GE P/N: 5269685)

Procedure Timing:

- 20 min for Old Style Filter
- 20 min for New Style Filter

Tools Required:

- Adjustable wrench
- Snap Ring Pliers (included in HEC or Cryogen Supply Header FRU Kit)

Old Style Filter Replacement:

1. Refer to the illustration below, and actuate the three yellow-handled valves in the order described:
 - a. Open lower right valve: Callout “A”
 - b. Close middle valve: Callout “B”
 - c. Close left valve: Callout “C”

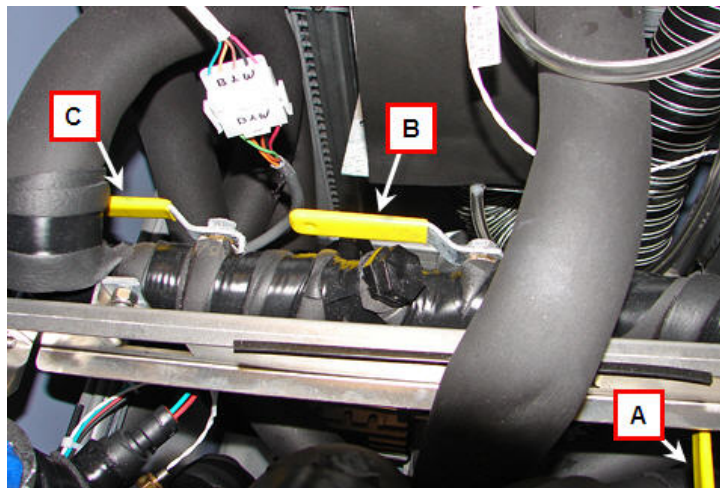


Illustration 1: Location of Cryocooler Compressor Valves

2. Remove the black plastic hex cap shown below.



Illustration 2: Location of Hex Cap over Filter

3. Pull out the wire mesh cylindrical screen from the exposed valve and replace it.
4. Proceed as follows to complete the process:
 - a. Install a clean filter and tighten down securely.
 - b. Replace the hex cap.
 - c. Open the left valve (C).
 - d. Open the middle valve (B).
 - e. Close the right valve (A).

New Style Filter Replacement:

1. Close the Jomar Filter Ball Valve in the header by turning the handle perpendicular to the pipe of the supply header. (Open and closed positions are shown below.)



Illustration 3: Filter Ball Valve Open



Illustration 4: Filter Ball Valve Closed

2. Unscrew the top of the ball valve using an adjustable wrench.

3. Remove the snap ring from the valve with the snap ring pliers to expose the strainer basket (The snap ring pliers are included in the HEC).

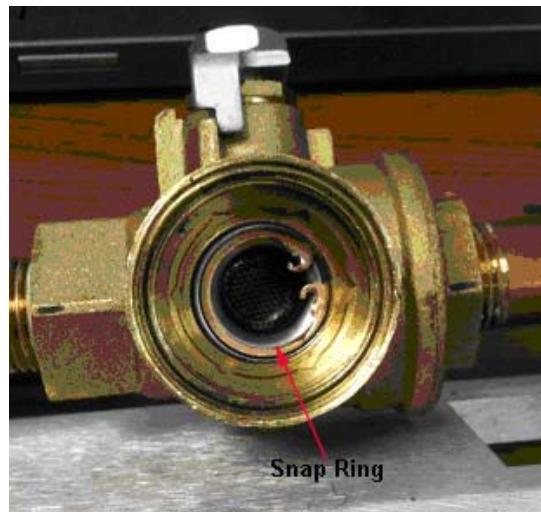


Illustration 5: Location of Snap Ring in Valve

4. Remove the strainer basket.



Illustration 6: Removing Strainer Basket from Valve

5. Replace the basket strainer and snap ring.
6. Replace the top of the ball valve, and tighten slightly past hand tight to ensure the O-ring is engaged.

REVISION HISTORY		
REV	ECN / DATE	REASON FOR CHANGE
A	6/11/2008	Initial Release
B	ECN# 17216 / 6/04/2010	Documenting changes requested by GE