

REVISION CONTROL INFORMATION

NOTE:

Before using this document please make sure that it is the latest revision using MyWorkshop.

This cover sheet controls the revision status of the entire document.

The revision of all sheets of this document will be at the same level.

Rev	Data	Author	Description
1	2011.Dec.05	Liu Wenzhang	Copy from 5371586 rev3
2	2013.Jun.25	Yang Jian	Add labels 5476711-xx and 5476712-xx
3	2013.Aug.28	Zhai Chengyao	Replace 5427868-x with 5454670-x

5443838TPH

SV 1.5T Multi-Language Label Installation Instructions

Technical Publication

Language Translation Key

Each label bears a base part number followed by a dash and a number that designates the language of that label. This key defines the language numerical suffixes on the labels. In this document the base label with a -x suffix will be shown. Replace x with the number for the desired language.

1: English	6: Italian	11: Portuguese	16: Czech	21: Finnish	26: Serbian
2: German	7: Dutch	12: Russian	17: Greek	22: Korean	27: Slovakian
3: French	8: Danish	13: Chinese	18: Bulgarian	23: Latvian	28: Turkish
4: Swedish	9: Norwegian	14: Hungarian	19: Croatian	24: Lithuanian	29: Slovenian
5: Spanish	10: Japanese	15: Polish	20: Estonian	25: Romanian	

1. PURPOSE

The purpose of this procedure is to increase the safety for the system operator and patient by applying the Safety Label. The safety labels, which must be applied at the appropriate location on the system, are available in different languages.

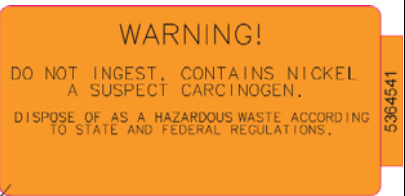
Service engineer must choose the safety label written in local language and apply it at the appropriate locations as described in this procedure.
If there is no space to put additional local language label, it is ok to put it on the English label with same meaning.

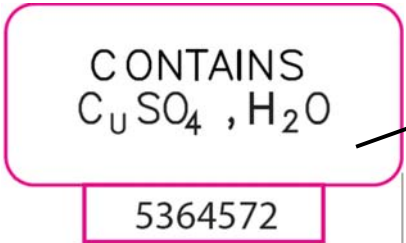
Note: The following English Label are already attached on product, there are no these Labels on English Label Kit.

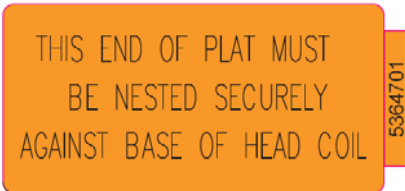
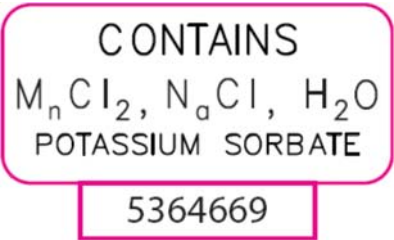
5370194-1	Laser Alignment Light Caution Label
5343503-1	LPCA cover warning label
5364539-1	RESET STEADY PULSE Label
5364540-1	Warning Label - contain nickel
5364541-1	Warning Label - contain nickel
5364572-1	Warning Label - contain CuSO4
5364669-1	Warning Label - contain MnCl2, NaCl
5364670-1	Warning Label - contain NiCl2, H2O
5364688-1	Warning Oxygen Levels Label
5364701-1	Nesting Plate Label
5366135-1	Host PC Drop Risk Warning Label
5366136-1	Authorized Personnel Only Warning Label
5366138-1	High Internal Temp Warning Label
5366141-1	P9180JU Caution Label

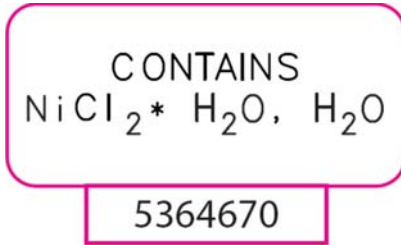
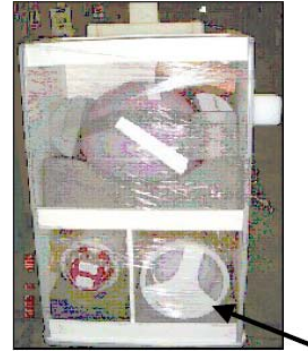
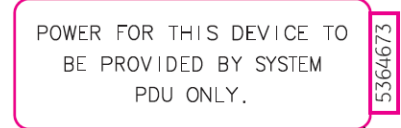
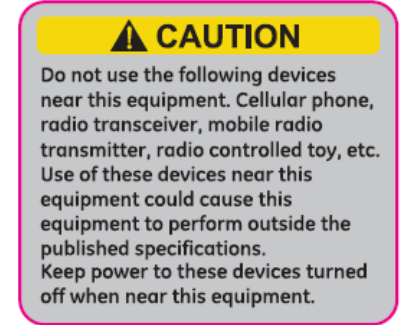
2. LABEL LIST

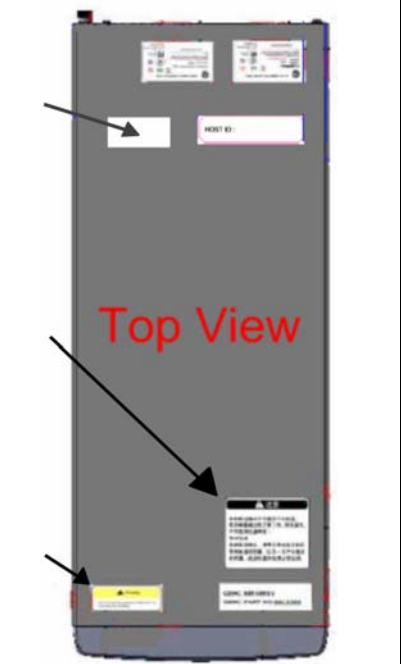
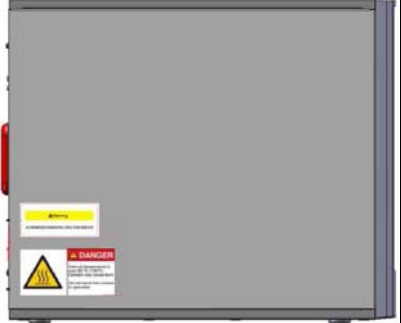
Description Content	Label Part Number	QTY	Label Image	Label Location on Part	Label Application Instruction	Location of part to be labeled
Label, Warning, Laser Radiation	5454670-X	1	 A yellow rectangular label with a black border. It features a black triangle with a sunburst symbol inside. The text reads: "LASER RADIATION", "DO NOT STARE INTO THE BEAM", "OR VIEW DIRECTLY WITH", "OPTICAL INSTRUMENTS", "CLASS 2M LASER PRODUCT", "1mW @ 655nm", "IEC 60825-1:2001", and "5454670" at the bottom left.	 A photograph showing the laser alignment light warning label applied to the front of a device, specifically over the magnet area.	Install the alternate language Laser Alignment Light Warning Label over the English labels.at the location shown	Front of magnet

Description Content	Label Part Number	QTY	Label Image	Label Location on Part	Label Application Instruction	Location of part to be labeled
MR Patient Alert – Patient Alert Labels(Black and white)	5364539-X	1			Open the desktop collector box and find the Patient Alert, Part #:46-317758P1. Install the appropriate 5364539-X alternate language label on the Patient Alert box as shown. Place the alternate language labels below the English labels. NOTE: The label should have RESET with the RED LED and STEADY/PULSE with the GREEN LED.	The Patient Alert is located inside the SV Desktop Collector.
MR Phantom Labels (orange and black)	5364540-X	2			Install the 5364540-X alternate language label on each 4.38" (100 mm) sphere. Place the new label beside the existing English label.	These small spheres are located in one of the following Mcats: M3335JA, M3335JB
MR DQA Phantoms - MR Phantom Labels (orange and black)	5364541-X	1			Install the alternate language 5364541-X label on the DQA/HL Phantom. Place the appropriate label on the end of the DQA phantom.	The DQA phantom is contained within the phantom cart. The phantom cart for 1.5T is 2133388-11. Remove the shrink-wrap and locate the DQA phantom on the phantom cart. Label it accordingly.

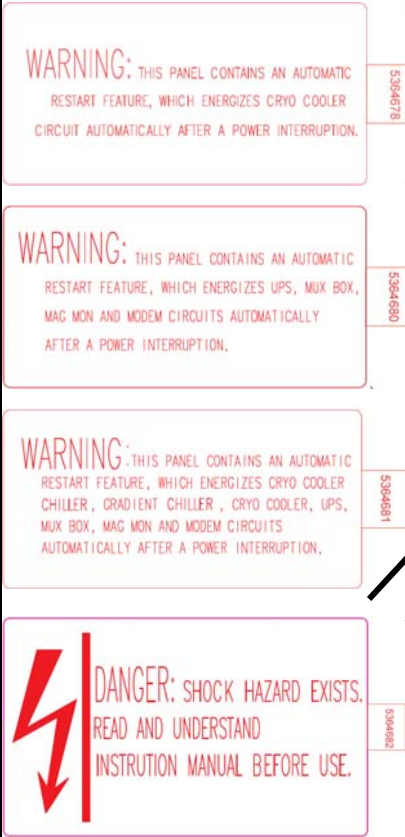
Description Content	Label Part Number	QTY	Label Image	Label Location on Part	Label Application Instruction	Location of part to be labeled
LV Shim phantoms - 1.5T only. MR Phantom Labels (clear and black)	5364572-X	3		  	Place the 5364572-X alternate language label on each piece of the LV Shim phantom beside the existing English labels. NOTE: Examine all the phantoms and the label sheet to verify that you have installed all the applicable labels.	Remove shrink-wrap and locate the LV Shim phantom assembly on the phantom cart and label accordingly. 

Description Content	Label Part Number	QTY	Label Image	Label Location on Part	Label Application Instruction	Location of part to be labeled
MR Nesting Plate in Phantom Cart – MR Phantom Labels (clear and black)	5364701-X	1		 5364580 5364701 5364700	Install the alternate language labels on the Short Nesting Plate beside the existing English label. Place the 5364580-X alternate language label on the Short Nesting Plate. Apply the alternate language 5364701-X label below the existing English label.	This Nesting Plate is mounted on the backside of the phantom cart. 
MR Body Loader-MR Phantom Labels (clear and black)	5364669-X	1			Place the 5364669-X alternate language label on the Body loader phantom beside the existing English label. NOTE: Examine all the phantoms and the label sheet to verify that you have installed all the applicable labels.	Remove shrink-wrap and locate the body / short loader phantom on the phantom cart and label accordingly. 

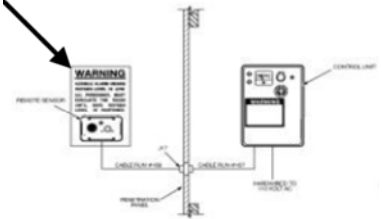
Description Content	Label Part Number	QTY	Label Image	Label Location on Part	Label Application Instruction	Location of part to be labeled
MR Short Loader Sphere – MR Phantom Labels (clear and black)	5364670-X	1			Place the 5364670-X alternate language label on the Short Loader sphere beside the existing English label. NOTE: Examine all the phantoms and the label sheet to verify that you have installed all the applicable labels.	Remove the shrink-wrap and locate the Short Loader sphere phantom on the phantom cart. Label it accordingly. 
MR Simplified OC - "Power for this device to be provided by System PDU only" Label (white and black)	5364673-X and 5366141-X and 5366135-X	1 each	 	Rear Side	Install the 5364673-X, 5366141-X, 5366135-X alternate language label beside the existing English labels on the lower rear of the Simplified OC.	MR simplified OC 5341380

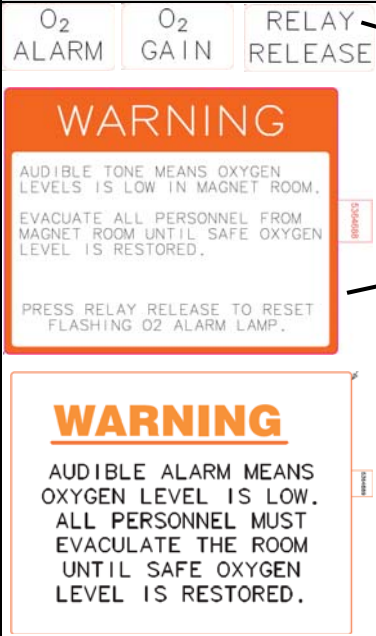
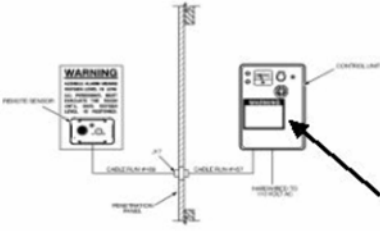
Description Content	Label Part Number	QTY	Label Image	Label Location on Part	Label Application Instruction	Location of part to be labeled
				 Top View Front Side		
Host PC High Temperature Warning Label	5366136-X and 5366138-X	1 each	 	 Left side OC chassis	Install the 5366136-X, 5366138-X alternate language label beside the existing English labels on the lower rear of the Simplified OC.	MR simplified OC 5341380

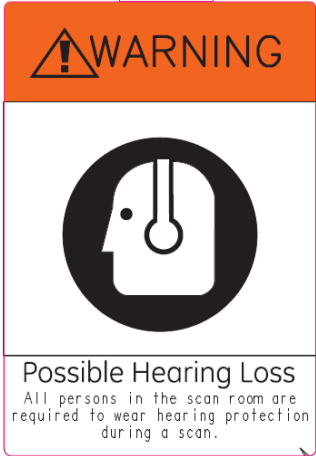
Description Content	Label Part Number	QTY	Label Image	Label Location on Part	Label Application Instruction	Location of part to be labeled
MR TRM MDP – MR MDP Labels (white and red)	5364676-X and 5364677-X and 5364674-X	1 each			Use the appropriate alternate language labels 5364676-X, 5364677-X, and 5364574-X for the MDP M3088TM. Install the labels beside the English labels. NOTE: The order may not include this Main Disconnect Panel (MDP).	The MDP M3088TM is in a large rectangular box. Cut the bands to open the box, and then apply the labels.

Description Content	Label Part Number	QTY	Label Image	Label Location on Part	Label Application Instruction	Location of part to be labeled
MR TRM MDP – MR MDP Labels (white and red)	5364678-X 5364680-X 5364681-X 5364682-X	1 each			Use the alternate language labels 5364678-X, 5364680-X, 5364681-X, 5364682-X for the MDP M3088TM. Install these labels beside the English labels. NOTE: The order may not include this Main Disconnect Panel (MDP)	The MDP M3088TM is a large rectangle box. Cut bands to open box, apply labels.

Description Content	Label Part Number	QTY	Label Image	Label Location on Part	Label Application Instruction	Location of part to be labeled
MR 1.5T MDP – MR 1.5T MDP Labels (white and red)	5364683-X and 5364684-X and 5364685-X	1 each			Install the alternate language labels 5364683-X, 5364684-X, 5364685-X on the MDP as shown. NOTE: The order may not include this Main Disconnect Panel (MDP).	The MDP Power Panel MR Signa Main, E4502SP is a stand-alone item.

Description Content	Label Part Number	QTY	Label Image	Label Location on Part	Label Application Instruction	Location of part to be labeled
Cryogen Cooler Compressor Warning Label (white and black)	5364686-X	1			Install the 5364686-X alternate language label below the existing English label.	 Cryogen Cooler Front and Rear
Remote Oxygen Alarm warning sign (orange, black, and white)	5364687-X	1			Mount the remote oxygen monitor with the warning sign on the RF shield wall within the scan room.	

Description Content	Label Part Number	QTY	Label Image	Label Location on Part	Label Application Instruction	Location of part to be labeled
Oxygen Monitor Labels (orange, black and white) Oxygen Monitor labels (clear and black)	5364688-X and 5364689-X and 5370458-X	1 each			Mount the alternate language oxygen monitor warning sign, 5364689-X, on the front of the Oxygen monitor. Mount the three small labels, 5364690-X, beside the existing English labels.	

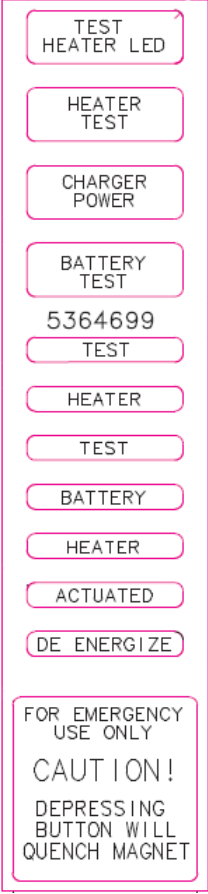
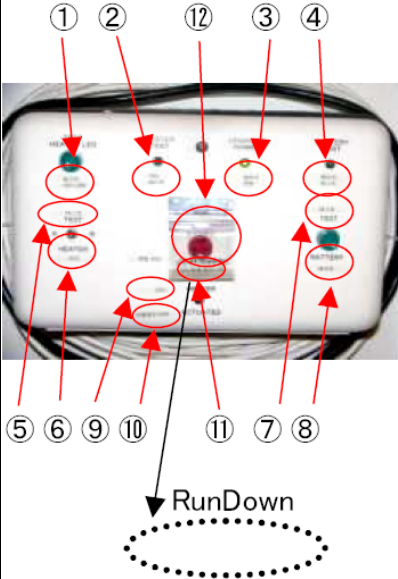
Description Content	Label Part Number	QTY	Label Image	Label Location on Part	Label Application Instruction	Location of part to be labeled
Hearing Loss Warning Label (orange, black, and white)	5364691-X	1	5364691 		Place the alternate language Hearing Loss Warning Label 5364691-X on the scan room door. Also described in the System Installation instructions.	See "Label Location on the Part" column.

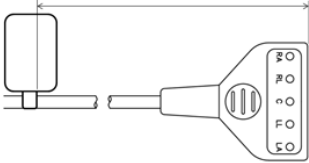
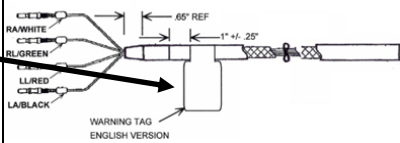
Description Content	Label Part Number	QTY	Label Image	Label Location on Part	Label Application Instruction	Location of part to be labeled
Authorized Personnel Warning sign (red, black, and white)	5364692-X	2			AUTHORIZED PERSONNEL WARNING. Set the sign at the entrance in front of the scan room door(s) when performing service procedures.	See "Label Location on the Part" column
Improper Coolant Warning Label (yellow and black)	5364693-X	1		 Water Heat Exchanger	Apply the alternate language Improper Coolant Warning Label 5364693-X on or near the Gradient/Coldhead Chiller Unit. The Chiller units are either M3088TL or M3088TK.	Determined by the site configuration.

Description Content	Label Part Number	QTY	Label Image	Label Location on Part	Label Application Instruction	Location of part to be labeled
Magnetic Warning Label 1.5" x 2.5" (orange, black, and white)	5364694-X	2		Outside of Scan Room	Place outside the scan room as needed.	Determined by the site configuration.
Magnetic Warning Label 3" x 5" (orange, black, and white)	5364695-X	10		Outside of Scan Room	Place outside the scan room as needed.	Determined by the site configuration.

Description Content	Label Part Number	QTY	Label Image	Label Location on Part	Label Application Instruction	Location of part to be labeled
5 Gauss Warning Sign (orange, black, and white)	5364696-X	10		Outside of 5 gauss line.	Place outside the 5 gauss line as needed.	Determined by the site configuration.
1.5T Magnet Room Door Warning Sign (orange, yellow, white and black)	5364697-X	2			Place the alternate language Magnetic Warning Sign 5364697-X on the scan room door. For 1.5T systems only.	See "Image of Label Location on the Part" Column.

Description Content	Label Part Number	QTY	Label Image	Label Location on Part	Label Application Instruction	Location of part to be labeled
Ferrous Components Warning (orange and black)	5364698-X	1			Install the alternate language label 5364698-X on the magnet dock assembly. Place the new label beside the existing English label. The Dock Assembly is located in the front of the magnet where the table attaches to the magnet.	Front of magnet

Description Content	Label Part Number	QTY	Label Image	Label Location on Part	Label Application Instruction	Location of part to be labeled
Magnet RunDown Unit Labels (clear and black)	5364699-X	1 each			Install the 5364699-X alternate language labels on the MRU. Place the new labels beside or below the English labels.	In Scan Room. Determined by site configuration.

Description Content	Label Part Number	QTY	Label Image	Label Location on Part	Label Application Instruction	Location of part to be labeled
MR Coil Connection on LPCA – Coil Connection Warning on Rear Ped or LPCA Label (orange and black)	5343503-X	1 each			Open the appropriate collector and install the applicable 5343503-X alternate language label onto the LPCA. Place the alternate language label beside the existing English label.	LPCA cover of HDsv Rear Pedestal 5330100 
MR Magnet ECG Port Cable Label	5476711-X*	1			If existing label has “Part No. E8819RD” DO NOT APPLY 5476711-X. Otherwise install the 5476711-X* label directly on the existing label.	Install according to the Image of “English Label Location on the Part” column
MR High Impedance ECG Cable Label	5476712-X*	1			If existing label has “Part No. 9218A” DO NOT APPLY 5476712-X. Install the 5476712-X* label directly on the existing label. Copy Date Code  and Lot Number  from English Label using a permanent marker.	Install according to the Image of “English Label Location on the Part” column