

REVISION CONTROL INFORMATION

NOTE:

Before using this document please make sure that it is the latest revision using MyWorkshop.

This cover sheet controls the revision status of the entire document.

The revision of all sheets of this document will be at the same level.

Rev	Date	Author	Description
1	2009.Nov.05	Jie Qing	Initially Created
2	2009.Dec.29	Jie Qing	Add a note on how to dealing with English label on field, if there have been English label on part; Add LPCA cover label Instruction
3	2010.Aug.31th	Jie Qing	• Add Slovenian language in table of “Language Translation Key” on preface • Change label quantity of 5364670-X to “1”, 5364693-X to “1”, 5364694-X to “2”, change 5364695-X and 5364696-X to “10”. Change PN of “5364598” to “5364698” because of typo. • Delete label 5364579, 5364580, 5364700 as there is no multi-language label for these three labels
4	2011.Aug.14th	Liu Wenzhang	Multiple changes made per ECR2115442
5	2013.Jun.25th	Yang Jian	Add labels 5476711-xx and 5476712-xx
6	2013.Aug.28	Zhai Chengyao	Replace 5427868-x with 5454670-x
7	2014.Sep.19	Zhai Chengyao	Replace 5116163-x with 5457590-x Replace 5116165-x with 5458376-x
8	2014.Nov.13	Zhai Chengyao	Add 5537470-9 Chinese laser warning label

5371586TPH

SV 1.5T Multi-Language Label Installation Instructions

Technical Publication

Notes:

1. Each sign/label has a base part number followed by a dash and a number that designates the sign/label language.

Examples:

5120663-2, French

5120663-16, Dutch

2. The “Label Image” column should be an image of the English base part number label.
3. The “Image of English Label Location on the Part” column may contain older revisions of the label.
4. In all cases, where there is room, the foreign language label shall be placed adjacent to the English language label.

1. PURPOSE

The purpose of this procedure is to increase the safety for the system operator and patient by applying the Safety Label.

The safety labels, which must be applied at the appropriate location on the system, are available in different languages.

Service engineer must choose the safety label written in local language and apply it at the appropriate locations as described in this procedure.

If there is no space to put additional local language label, it is ok to put it on the English label with same meaning.

2. LABEL LIST

Description Content	Label Part Number	QTY	Label Image	Label Location on Part	Label Application Instruction	Location of part to be labeled
Label, Warning, Laser Radiation	5454670-x	1	 A yellow rectangular label with a black border. It features a black triangle with a sunburst symbol inside. The text reads: "LASER RADIATION", "DO NOT STARE INTO THE BEAM OR VIEW DIRECTLY WITH OPTICAL INSTRUMENTS", "CLASS 2M LASER PRODUCT", "1mW @ 655nm", "IEC 60825-1:2001", and the part number "5454670" at the bottom left.	 A photograph showing the English warning label applied to the front cover of the magnet bore of an MRI machine.	Install the alternate language Laser Alignment Light Warning Label over the English labels at the location shown (except for product shipped to China).	On front cover of magnet.  A photograph showing the interior of the MRI magnet bore, with the warning label visible on the front cover.
LASER ALIGNMENT LIGHT ED 2.0 WARNING LABEL - Chinese	5537470-9	1	 A yellow rectangular label with a black border. It features a black triangle with a sunburst symbol inside. The text is in Chinese: "激光辐射", "勿直视或通过光学仪器直接观看光束", "2M类激光产品", "波长650 ± 10nm", "最大输出功率小于1毫瓦", "GB7247.1-2012", and the part number "5537470-9" at the bottom left.	 A photograph showing the Chinese warning label applied to the front cover of the magnet bore of an MRI machine.	Install the Chinese language Laser Alignment Light Warning Label over the English labels at the location shown for product shipped to China.	On front cover of magnet.  A photograph showing the interior of the MRI magnet bore, with the warning label visible on the front cover.

Description Content	Label Part Number	QTY	Label Image	Label Location on Part	Label Application Instruction	Location of part to be labeled
MR Patient Alert – Patient Alert Labels	5224905-x	1			Find the Patient Alert, Part # 46-317758P1, on the operator workspace area. Install the appropriate 5224905-x alternate language label on the Patient Alert box as shown. Place the alternate language labels below the English labels. NOTE: The label should have RESET with the RED LED and STEADY/PULSE with the GREEN LED.	The Patient Alert is on the operator workspace area. 
MR Phantom Labels	5229689-x	2			Install the alternate language warning label on each 4.38" (100 mm) sphere. Place the new label adjacent to the existing English label..	These small spheres are located with the customer MR system calibration phantoms.
MR DQA Phantoms - MR Phantom Labels	5229690-x	1			Install the alternate language 5229690-x warning label on the DQA/HL Phantom. Place the appropriate label on the end of the DQA phantom.	The DQA phantom is contained within the optional phantom cart (shown below) or is located within the customer MR system phantom storage area.

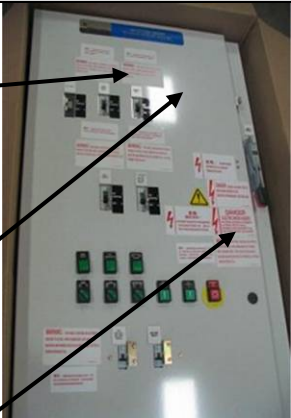
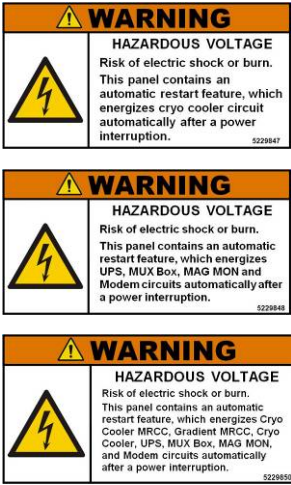
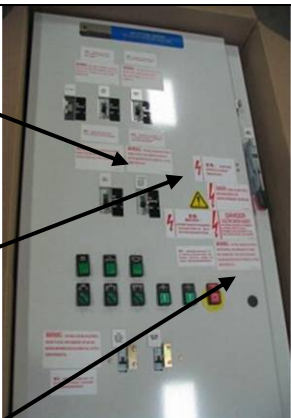


Description Content	Label Part Number	QTY	Label Image	Label Location on Part	Label Application Instruction	Location of part to be labeled
						
LV Shim phantoms - 1.5T only. MR Phantom Labels	5229702-x	3	 WARNING DO NOT INGEST Contains CuSO₄, H₂O. Dispose of as a hazardous waste according to local regulations. See GE MSDS # WM-61. 5229702		Install the alternate language warning label on each piece of the LV Shim phantom adjacent to the existing English labels. Note: This Phantom is only used on 1.5T Systems. NOTE: Examine all the phantoms and the label sheet to verify that you have installed all the applicable labels.	The LV Shim phantom assembly is located on the optional phantom cart (shown below) or is located within the customer MR system phantom storage area. 

Description Content	Label Part Number	QTY	Label Image	Label Location on Part	Label Application Instruction	Location of part to be labeled
MR Nesting Plate in Phantom Cart – MR Phantom Labels	5229691-x	1		 5229691-x	Install the alternate language labels on the Short Nesting Plate adjacent to the existing English label. Apply the alternate language 5229691-x label below the existing English label.	This Nesting Plate is mounted on the backside of the phantom cart.  
MR Body Loader-MR Phantom Labels	5229703-x	1			Install the alternate language warning label on the Body loader phantom adjacent to the existing English label. NOTE: Examine all the phantoms and the label sheet to verify that you have installed all the applicable labels.	The body / short loader phantom is located on the optional phantom cart (shown below) or is located within the customer MR system phantom storage area..

Description Content	Label Part Number	QTY	Label Image	Label Location on Part	Label Application Instruction	Location of part to be labeled
						
MR Short Loader Sphere – MR Phantom Labels	5229689-x	2			Install the alternate language warning label on the Short Loader sphere adjacent to the existing English label. NOTE: Examine all the phantoms and the label sheet to verify that you have installed all the applicable labels.	The Short Loader sphere phantom is located on the optional phantom cart (shown below) or located within the customer MR system phantom storage area. 

Description Content	Label Part Number	QTY	Label Image	Label Location on Part	Label Application Instruction	Location of part to be labeled
MR Simplified OC - "Power for this device to be provided by System PDU only" Label	5427603-X and 5427022-X	1 each			Install the 5427603-X, 5427022-X alternate language label on the top of the Simplified OC.	MR simplified OC 5341380
Host PC High Temperature Warning Label	5427020-X	1			Install the 5427020-X alternate language label on the top of the Simplified OC.	MR simplified OC 5341380

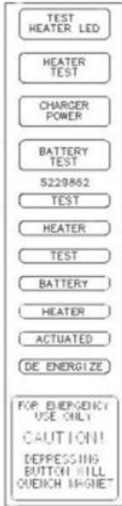
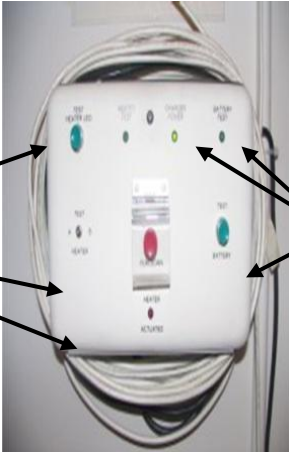
Description Content	Label Part Number	QTY	Label Image	Label Location on Part	Label Application Instruction	Location of part to be labeled
MR TRM MDP (Main Disconnect Panel) - MR MDP Labels	5229844-x AND 5229845-x AND 5229846-x	1 each			Install the appropriate alternate language safety labels on the MDP M3088TM, adjacent to the English labels. Applicable to 1.5T and 3.0T TRM only. NOTE: The order may not include this MDP (Main Disconnect Panel).	The MDP is usually mounted on the wall in the equipment room as shown to the left.
MR TRM MDP (Main Disconnect Panel) - MR MDP Labels	5229847-x AND 5229848-x AND 5229850-x	1 each			Install the alternate language warning labels on the MDP adjacent to the English labels. Applicable to 1.5T and 3.0T TRM only. NOTE: The order may not include this MDP (Main Disconnect Panel).	The MDP is usually mounted on the wall in the equipment room as shown to the left.

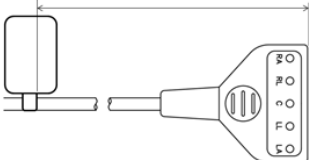
Description Content	Label Part Number	QTY	Label Image	Label Location on Part	Label Application Instruction	Location of part to be labeled
Cryogen Cooler Compressor Warning Label	5229863-x	1			Install the alternate language warning label below the existing English label.	Cryogen Cooler Front and Rear
Remote Oxygen Alarm Warning Sign	5116164-x	1			Install the remote oxygen monitor with the warning sign on the wall within the scan room.	
Oxygen Monitor Label	5231683-x	1			Install the alternate language warning sign on the front of the oxygen monitor in the operator room.	

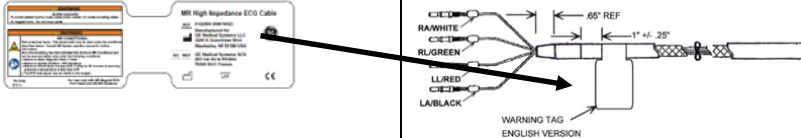
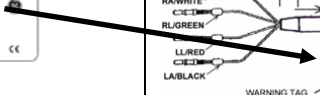
Description Content	Label Part Number	QTY	Label Image	Label Location on Part	Label Application Instruction	Location of part to be labeled
Authorized Personnel Warning Sign	5116158-x	2	 The sign features a yellow header with a black exclamation mark and the word "WARNING". Below this is a red circle with a diagonal slash over a black silhouette of a hand. At the bottom, the text reads: "HAZARDOUS AREA Authorized personnel only beyond this point. Equipment installation or servicing in progress."	N/A	Set the warning sign at the entrance in front of the scan room door(s) when performing service procedures.	See "Image of Label Location on the Part" Column
Improper Coolant Warning Label	5116159-x	2	 The label has a yellow header with the word "CAUTION" in black. Below it, the text states: "Improper coolant can cause serious damage to equipment." At the bottom, there is a blue circular icon of a person reading a manual, with the text "Refer to service manual for requirements" next to it.	 A photograph of a large, industrial-grade chiller unit with a metal grille on the left and a control panel on the right.	Install the alternate language caution label, on or near the Gradient/Coldhead Chiller Unit.	Determined by the site configuration.
Magnetic Warning Label 1.5" x 2.5"	5116161-x	10	 The label has a yellow header with a black exclamation mark and the word "WARNING". Below this is a red circle with a diagonal slash over a black silhouette of a hand holding a magnet. At the bottom, the text reads: "This equipment must not be taken into the magnet room. Damage to the equipment, MR system and personal injury could result."	Outside of Scan Room	Outside the scan room, install warning labels as needed on ferrous equipment that may accidentally enter the scan room such as chairs, hospital equipment, janitorial equipment, etc.	Determined by the site configuration.

Description Content	Label Part Number	QTY	Label Image	Label Location on Part	Label Application Instruction	Location of part to be labeled
5 Gauss Warning Sign	5457590-x	10		Outside of 5 gauss line.	Install warning signs outside the 5 gauss line as needed. Recommend posting these signs two days prior to the activation of the magnet.	Determined by the site configuration.
1.5T Magnet Room Door Warning Sign	5458376-x	2			Install the 1.5T Magnetic Warning Sign on the scan room door. If door opens outward, place signs on both sides of door. For 1.5T systems only.	See "Image of English Label Location on the Part" Column.
errous Components Warning	5327239-x or 5229860-x	1			Install the alternate language warning label on the magnet dock assembly. Place the new label adjacent to the existing English label. The Dock Assembly is located in the front of the magnet where the table attaches to the magnet.	At front of magnet on the floor.



Description Content	Label Part Number	QTY	Label Image	Label Location on Part	Label Application Instruction	Location of part to be labeled
						
Magnet Rundown Unit Labels	5229862-x	1 each			Install the alternate language labels on the MRU. Place the new labels adjacent to or below the English labels.	In Scan Room. Determined by site configuration.

Description Content	Label Part Number	QTY	Label Image	Label Location on Part	Label Application Instruction	Location of part to be labeled
MR Patient Handling for BRM only – MR Patient Handling Label	5255438-x	1		Blower shown with the cover off. 	Install the alternate language warning label, part number, on the top of the Blower Cabinet. Note: This is applicable to the 1.5T BRM only.	BRM Blower is mounted, usually on the floor, in the scan room outside the 200 gauss line with the cover on as shown below: 
MR Coil Connection on LPCA – Burn Hazard	2146831-x	1			Open the appropriate collector and install the applicable 2146831-x alternate language warning label onto the LPCA. Place the alternate language label adjacent to the existing English label.	Rear Pedestal LPCA Rear Bridge Assembly 
MR Magnet ECG Port Cable Label	5476711-x*	1			If existing label has "Part No. E8819RD" DO NOT APPLY 5476711-x. Otherwise install the 5476711-x*	Install according to the Image of "English Label Location on the Part" column

Description Content	Label Part Number	QTY	Label Image	Label Location on Part	Label Application Instruction	Location of part to be labeled
					label directly on the existing label.	
MR High Impedance ECG Cable Label	5476712-x*	1	 The label image shows two parts: an English label on the left with regulatory information and a CE mark, and a diagram on the right. The diagram illustrates the cable's construction with four conductors labeled RA/WHITE, C/WHITE, RL/GREEN, and LL/RED. A black conductor is labeled LA/BLACK. Dimensions shown include a 0.05" REF for the conductor spacing and a 1" +/- .25" for the overall width. A warning tag labeled 'WARNING TAG ENGLISH VERSION' is attached to the cable.	 A black arrow points from the English label in the 'Label Image' column to the 'WARNING TAG ENGLISH VERSION' on the cable diagram in the 'Label Location on Part' column, indicating where to apply the new label.	If existing label has "Part No. 9218A" DO NOT APPLY 5476712-x. Install the 5476712-x* label directly on the existing label. Copy Date Code  and Lot Number  from English Label using a permanent marker.	Install according to the Image of "English Label Location on the Part" column