

POWER UP / DOWN PROCEDURE

This section describes the system start-up/shutdown procedures. **Power-UP** is accomplished via the keyswitch on the Display Tower. **Power DOWN** of the entire system is accomplished via the General Utilities GUI. This is accessed via the “Service” option from the opening bootup GUI. This procedure covers the Service Application mode. *Scanning Application mode procedure is in the OP manual.*

SYSTEM POWER UP

1. Make sure the Gantry keyswitch is in the normal position. Turn the Display Tower keyswitch to the **ON** position. After the system boots, you will be prompted with the startup application display (or “4 Button Screen.”) You have **4 boot options**. Refer to Figure 62.:
A.) Select Clinical Scanning (or do nothing) and the system will start and connect to the acquisition server. The acquisition server will then power up the AcQImage Cabinet, and if that is done successfully, the acquisition server will power up the gantry and X-ray system.
B.) Select the Service Application option (within 15 seconds) by clicking on the Service button with the mouse pointer.
C.) Reboot the system. **D.)** Start the Diagnostic GUI (within 15 seconds) by clicking once in the bottom left corner, then entering PICKER <CR>. (All caps). The AcQImage board LEDs will display a power up test sequence. Refer to the [LED Interpretations](#).

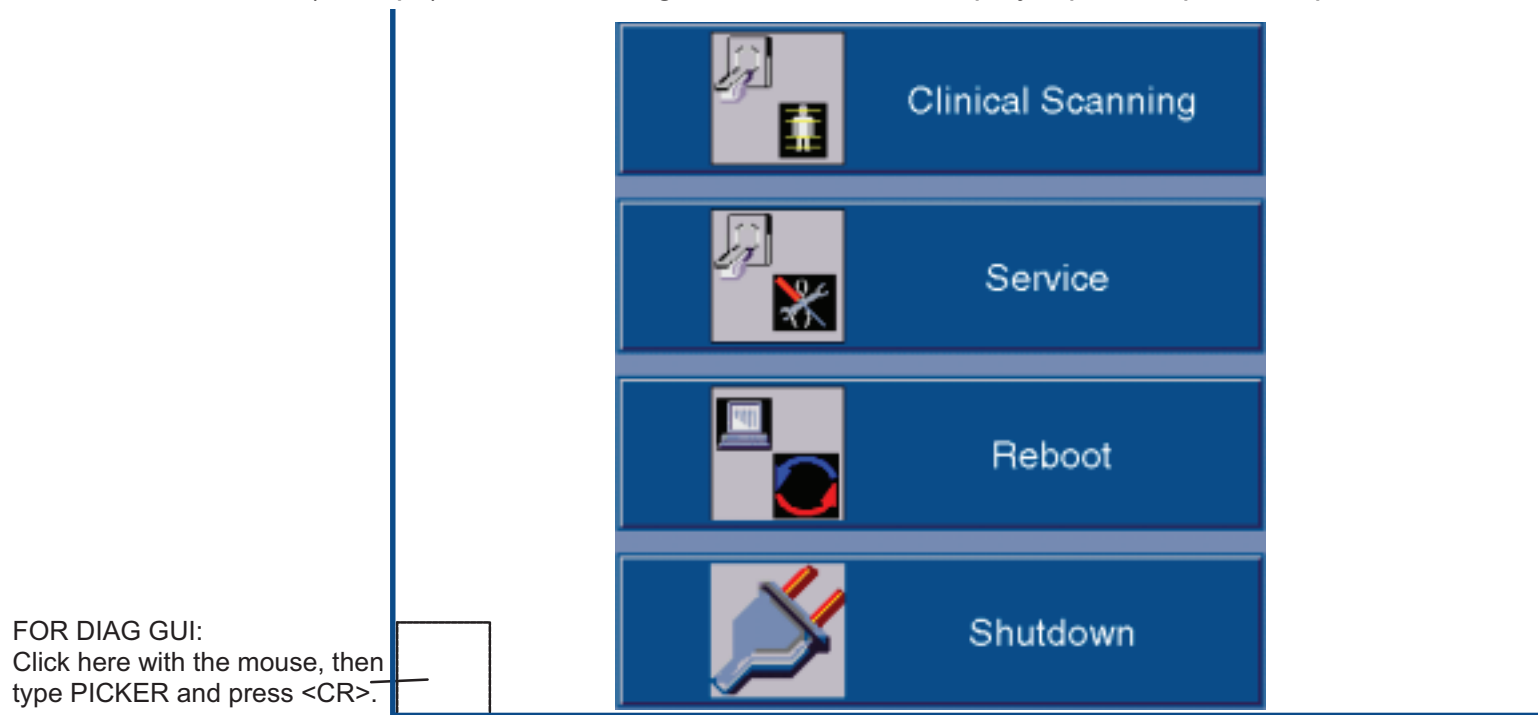


FIGURE 62 4 BUTTON SCREEN

PREPARE THE SYSTEM FOR USE (X-Ray Tube Warmup)

2. The X-Ray system needs to be properly warmed up prior to executing a scan procedure. X-Ray System Warmup must be completed at the start of every scan day or when the tube heat is below 6%. Tube heat is displayed in the Tube Warmup window.
3. To perform a tube warm-up, click on the Tube Warm Up button with the mouse pointer. Refer to Figure 62. This will initiate the warm-up sequence and **X-RAYS WILL BE EMITTED**. The tube warm-up status is displayed in the Tube Warmup Window. Refer to Figure 63.

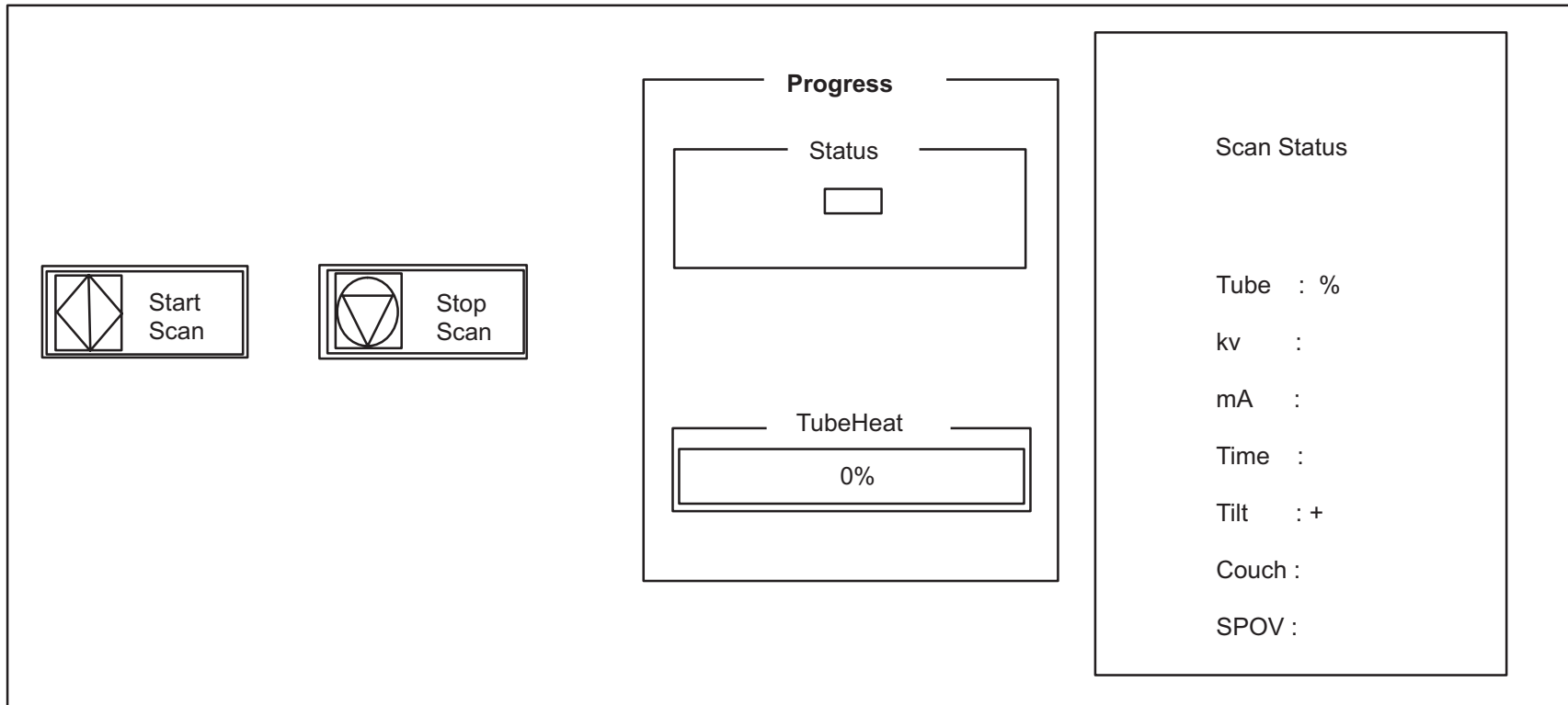


FIGURE 63 TUBE WARM-UP STATUS WINDOW

SYSTEM POWER DOWN

1. From the 4 Button Screen (Figure 62) click on SHUTDOWN with the mouse pointer.



CAUTION

DO NOT POWER DOWN THE SYSTEM WITH THE TUBE HEAT AT OR ABOVE 35%. FAILURE TO COMPLY MAY RESULT IN SHORTENED TUBE LIFE AND/OR ERRATIC TUBE BEHAVIOR.

System shutdown ensures that all images on disk are properly closed, all image maintenance functions properly terminated, the gantry parked, the X-Ray system sufficiently cooled and all films printed.

There are three shutdown selections: Shutdown Gantry, Shutdown Console, Shutdown System and Switch to Console as seen in Figure 64.

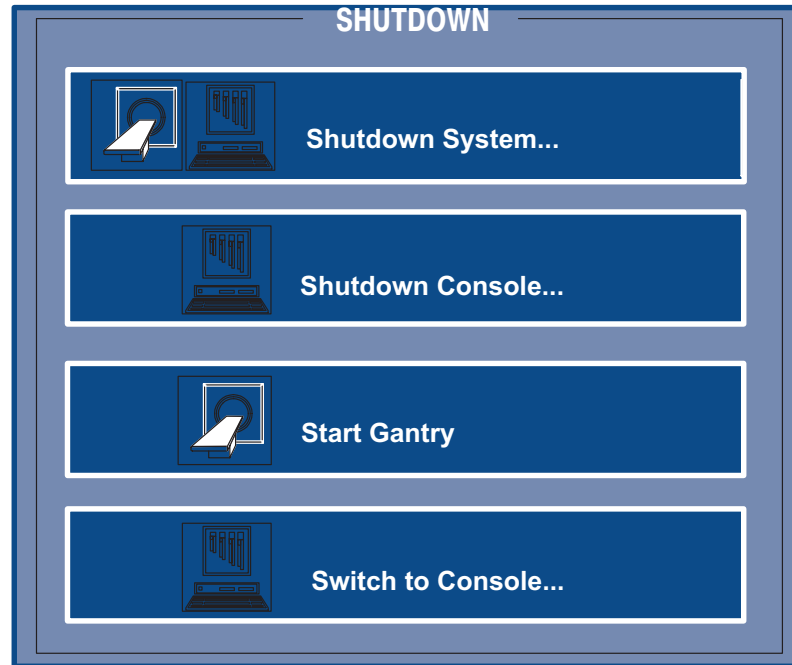


FIGURE 64 SHUTDOWN SCREEN

Shutdown Gantry: This option should be used, if possible, in place of the e-stop switches. An orderly shutdown of the gantry will be performed. Display Tower and/or AcQImage Cabinet operation may continue after selected. To cancel shutdown, click the Abort button.

- Shutdown System: Performs an orderly shutdown of the entire scanner system. To cancel shutdown, click the Abort button.
- Shutdown Console: This option shuts down only the console units. An orderly shutdown of the image archiver, filming and file system will be done.
- Start Gantry: Toggles Start/Shutdown Gantry depending on status of Gantry power.
- Switch to Console: Changes the GUI to the Console Application (Operator) screen.



WARNING



THE GUI SHUTDOWN REMOVES POWER FROM THE SYSTEM, BUT NOT AT THE INCOMING POWER LINES. DANGEROUS VOLTAGES ARE EXPOSED WHEN THE GANTRY COVERS ARE REMOVED. BEFORE SERVICING EQUIPMENT, REMOVE POWER FROM THE WALL BOX. FAILURE TO COMPLY CAN CAUSE SERIOUS INJURY OR DEATH TO PERSONNEL.

2. When prompted via a dialog box, turn the Display Tower keyswitch to the **OFF** position.
3. On the AcQImage Cabinet and Display Tower, ensure the rear breaker switch is turned off.

AcQImage Boards – Power Up Test LED Interpretations

AP, PB, GAM Board Diagnostic Flash Version 1.1.2

0 Hardware Initialization
1 CPU Version Check
2 Diagnostic Flash CRC Check
3 Operating System Flash CRC Check
4 CPU Memory Test
5 Data Cache Test
6 Unexpected Exception Test
7 IBC Test
8 NVRAM Single Location Test
9 PCI Bus Select Test
A Ethernet Loopback Test
B UART Loopback Test
C apmem for AP, pbmem for PB, gamem for GAM

C CIO Test (CP Board Diagnostic Flash Version 1.1.2)
D Synchronous Serial Chip Loopback Test (CP Board Diagnostic Flash Version 1.1.2)