



GE Medical Systems

Technical Publications

Direction 2389159

Revision 8

Signa[®] EXCITE 3.0T Field Replaceable Units (FRU)

Copyright © 2003 – 2009 by General Electric Company, Inc.

All Rights Reserved



GE Medical Systems

GE Medical Systems; Telex 3797371
P.O. Box 414, Milwaukee, Wisconsin 53201 U.S.A.
(Asia, Pacific, Latin America, North America)

GE Medical Systems – Europe; Telex 261794
Shortlands, Hammersmith, London W6 8BX U.K.



**COVER
BOTTOM
CROP LINE**

**COVER
BOTTOM
CROP LINE**

**5/16" RADIUS
CORNER TRIM
TYPICAL**



**COVER
LEFT
CROP LINE**

**COVER
FOLD
LINE**

**(2) CENTER STAPLES
FOR BINDING**

**COVER
RIGHT
CROP LINE**



**5/16" RADIUS
CORNER TRIM
TYPICAL**

REVISION HISTORY

REV/DATE REASON FOR CHANGE

0	Oct 1, 03	Initial release
1	Mar 8, 04	Added front end bell and LPCA parts. Added AirSys water chillers. Updated EXCITE software info, ACGD/PDU parts System cabinet parts & modules. Added fuses to RF Amp cabinet, modules to TAC, parts to phantom cart, NV Array Coil, and Option keys.
2	Nov 10, 04	Added cable replacement information. Revised Service Tools and Options. Added GOC AA to OW. Added BrainWave parts. Revised other part information.
3	Dec 30, 05	Added TRM and SRI parts. Updated OW monitor info, added LAN Hub part #. Updated RRF module part #'s. Changed NB RF Amp Cabinet FRU. Added alternate TAC vacuum pump module. Updated LV Shim phantom # and added DQA location bracket #.
4	Mar 13, 07	Added note for referring to 2138247 for parts not covered in this manual.
5	Jun 15, 07	Added LPCA LED harness. Updated ACGD SGA power supply from 2272460-2 to 2317092-2.
6	May 28, 08	Added RF Shield Change Out Kit. Rear Pedestal: updated Longitudinal Encoder part number (PQR 13175398) and added Encoder board part # 46-226552G2. MGD: Updated SCP2 Bd part #. Added 5181495-3 MGD P.S. Coils: Added T/R Extremity Coil (PQR 13102505) and updated NV Array Coil part #'s (PQR 13113173).
7	Jun 15, 09	Added Flexible Coupling (PQR 13206931). Updated Monitor offerings on OW (PCN 290188). Removed DVD-RW kit from OW (PQR 13247342). Also added 73Gb HDD (ECR 2084635). In RRF Chassis, added 5250028-2 RF-DIF (ECO 2080729).
8	Nov 3, 09	Added new GOC AA 2395043-2 (ECO 2084000)

SEND COMMENTS ABOUT THIS MANUAL TO:

GE MEDICAL SYSTEMS
MANAGER, SERVICE ENGINEERING
W-808
3200 N. GRANDVIEW BLVD.
WAUKESHA, WISCONSIN 53188
USA

This manual belongs to:

INTRODUCTION TO FIELD REPLACEABLE UNIT (FRU)

This direction is a guide to the FRU's.

Parts designated as a FRU are critical to maintain the safety and/or functionality of the system. FRU's are available on a Priority 1 basis from Service Parts Operation.

As a reference aid, the FRU's are listed under an upper-level part number where they are used. In most cases the upper-level part number is not a FRU. In some cases the parts listed under an upper-level assembly may have been nested in intermediate sub-assemblies which are not shown.

Where possible, FRU's which are part of Option Catalogs are shown under the area where they would be installed, rather than being listed in the Coils or Options Sections (SEE TABLE OF CONTENTS).

LIST OF EFFECTIVE PAGES


<u>PAGE</u>	<u>REV</u>	<u>PAGE</u>	<u>REV</u>	<u>PAGE</u>	<u>REV</u>
Cover Page ..	8	10 to 11	0	22	6
1	8	12 to 13	1	23	7
2	0	14	0	24	0
3 to 4	3	15	7	25 to 26	2
5	6	16	8	27	3
6	3	17 to 18	0	28	1
7	7	19	5	29 to 33	3
8	5	20	0	34	6
9	4	21	1	35 to 36	2

2


NOTES

TITLE	PAGE
INTRODUCTION TO FIELD REPLACEABLE UNIT(FRU)	1
LIST OF EFFECTIVE PAGES	1
TABLE OF CONTENTS	3
PICTORIAL TABLE OF CONTENTS	4
MAGNET AND ENCLOSURE (ASSEMBLIES)	5, 6, 7, 8, & 9
MISCELLANEOUS LIST (OF PARTS NOT SHOWN)	5
BRM & TRM LIST (OF PARTS NOT SHOWN)	5
FRONT COSMETIC COVER & PARTS (BACK VIEW)	5
ENCLOSURE & VACUUM RELATED	6
FRONT & REAR ENDBELL	6
COIL ASSEMBLY & BRIDGE SUPPORT	6
LONGITUDINAL MOTOR DRIVE	7
REAR PEDESTAL PARTS	7
COIL ADAPTERS	7
LOW PROFILE CARRIAGE ASSEMBLY	8
BRIDGE	8
TOUCH-UP PAINTS	9
PATIENT SUPPORT TABLE	10
CHILLER	12
PENETRATION PANEL	13
OPERATOR WORKSPACE – EXCITE	15
OPERATOR WORKSPACE COMPUTER CABINET – LINUX	16
WORKSPACE SCIM AND KEYBOARDS	17
MAIN DISCONNECT PANEL	18
ACGD/PDU CABINET	19
SYSTEM CONTROL CABINET	21, 22 & 23
RF ACCESSORY CABINET	25 & 26
NB RF AMP CABINET	27
8KW MNS BB RF AMP CABINET (OPTIONAL)	28
TWIN ACCESSORY CABINET (TAC)	29 & 30
REPLACEMENT CABLES	31
SERVICE TOOLS	32
ACCESSORIES	34
COILS	34
EXCITE OPTION KEYS	35
BrainWave OPTIONS	36

INTRO TO FRU'S	PAGE 1
CHILLER PARTS	PAGE 12
CABLES	PAGE 31
ACCESSORIES	PAGE 34
COILS	PAGE 34
OPTION KEYS	PAGE 35
BrainWave	PAGE 36

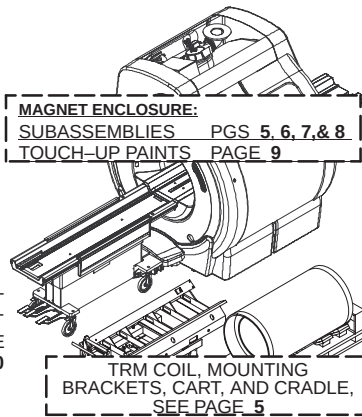

2241999
 PATIENT AIR SUPPLY
 BLOWER CABINET (MG6)

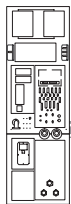


PENETRATION
 PANEL
 PAGE 13

PATIENT
 SUPPORT
 TABLE
 PAGE 10

MAGNET ROOM

OPERATOR/EQUIPMENT ROOM


EXCITE II OPERATOR
 WORKSPACE WITH LINUX
 PAGES 15, 16, & 17

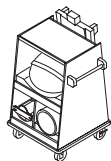

2219400
 MAGNET
 MONITOR
 (MSM)




MAIN
 DISCONNECT
 PANEL
 PAGE 18

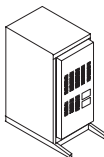


ACCESSORY
 CABINET
 PAGES 29 & 30




SERVICE TOOLS
 PAGE 32

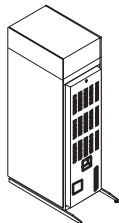


ACGD
 CABINET
 PAGE 19

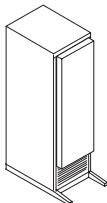


SYSTEM
 CABINET
 PAGES
 21, 22, & 23

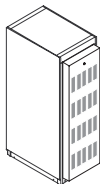


RF ACCESSORY
 CABINET
 PAGES
 25 & 26

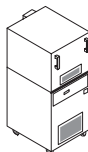

NB RF AMP
 CABINET
 PAGE 27

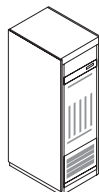

8KW BB RF
 AMP CABINET
 PAGE 28

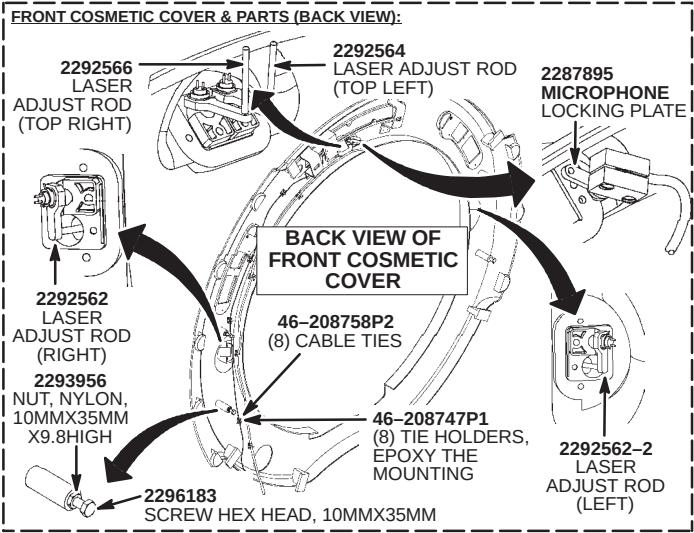


NOTE: The following LCC 3.0T FIXED W/TRM & Enclosure related parts Not Shown;
MISCELLANEOUS:


2203803	SPLIT O-RING
2242756	RELIEF VALVE, 3/8 NPT X 7 PSIG, GEN
2299843	MAGNET PRESSURE TRANSDUCER, SUNX DP2-41N
46-294231G1	MAGNET RUNDOWN UNIT
46-294231P1	BUSSMAN, F1, 110VAC .400 GDC
46-264758G1	REMOTE INTERCOM BOARD, MG3 A2 A1
TRM:	
2144093-9	TRM WITH MOUNTING BRACKETS, CART, AND CRADLE
2218148	COLD SHRINK (FOR REAR ENDBELL AND TRM COIL)
2208999-7	RF BODY COIL AND CASE
2299498	COOLING CIRCUIT REPAIR KIT
5113878	TRM RF SHIELD
2389049	RF SHIELD CHANGE OUT KIT
2293612	MALE CONNECTOR FOR 3/4" HOSE TO 3/4" NPT
5111938	BORE TEMPERATURE SENSOR
2383070	DIODE SLOT BOARD
or	
5250104	DIODE SLOT BOARD

ENCLOSURE & VACUUM RELATED PARTS:

2293771
RADIAL SHIM SPACER, 1.5MM
2293771-2
RADIAL SHIM SPACER, 0.8MM
2241176
RUBBER ISOLATION PAD

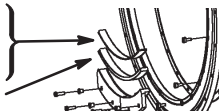**FRONT ENDBELL PARTS:
(TRM ONLY)**


2295639
O-RING,
ACCESS
PANELS

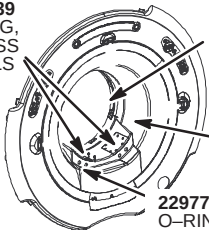

GASKET,
ENDBELL
TO RF TUBE
2294358, 4mm
2294358-2, 6mm
2294358-3, 9.5mm
2294358-4, 8mm

FRONT
ENDBELL ASM
2296177

2297761
O-RING, WATERLINE
FEEDTHROUGH

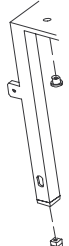
**BRIDGE SUPPORT PARTS:
(TRM ONLY)**

2293539
BRIDGE
SUPPORT
BUSHING

2293956
NYLON HEX NUT,
M10 X 1.5MM
X 9.5MM HIGH

**REAR ENDBELL ASSEMBLY
(TRM ONLY)**

2284170
REAR ENDBELL
VACUUM ASSEMBLY;
& 2218148
COLD SHRINK
(FOR REAR ENDBELL &
TRM COIL)
& 2218149
COMPRESSION
SPLICER
(FOR REAR ENDBELL)

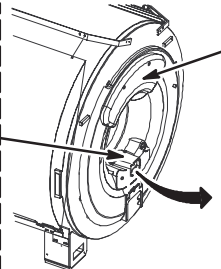

2358470
COMFORT MODULE
ASSEMBLY

**GRADIENT GASKETS:
(TRM ONLY)**




2297540
(4MM)
2297540-2
(6MM)

NOTE: The following parts are Not Shown.


2354050	SPECTRO T/R ASM (OPTION)
2381249	MNS LPCA CABLE HARNESS (OPTION)
2297760	VACUUM GREASE, MSDS# 8366442
2297672	GLYCOL, 1 GALLON/3.8 LITERS



NOTE: The following REAR PEDESTAL parts Not Shown:

2362995-2	8 CHANNEL SWITCH ASM
5115910	SRI 2 MODULE
or	
5121793	SRI 3 UPGRADE KIT

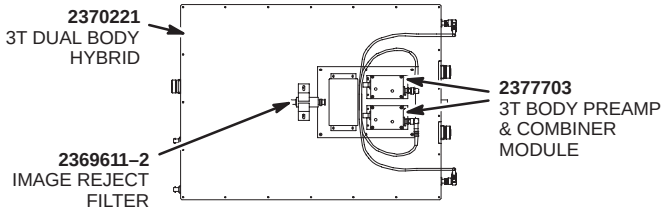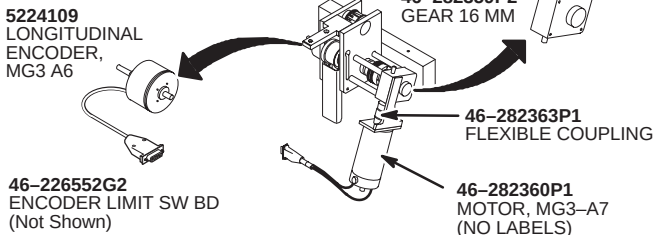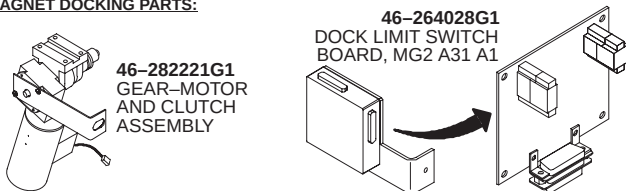3.0T HYBRID SPLITTER:**LONGITUDINAL MOTOR DRIVE PARTS:****MAGNET DOCKING PARTS:****REMOTE INTERCOM BOARD**

6250012
REMOTE INTERCOM
BOARD, MG3 A2 A1

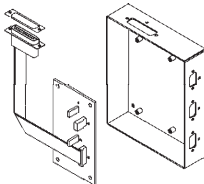


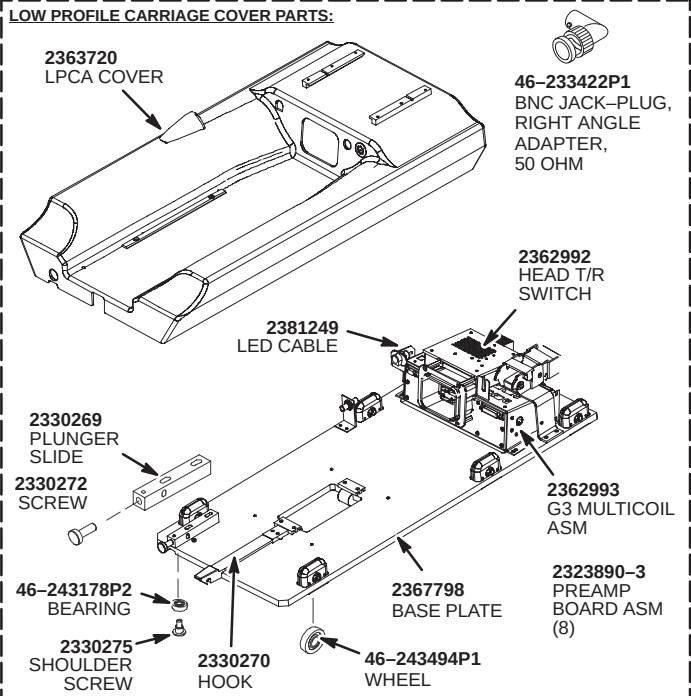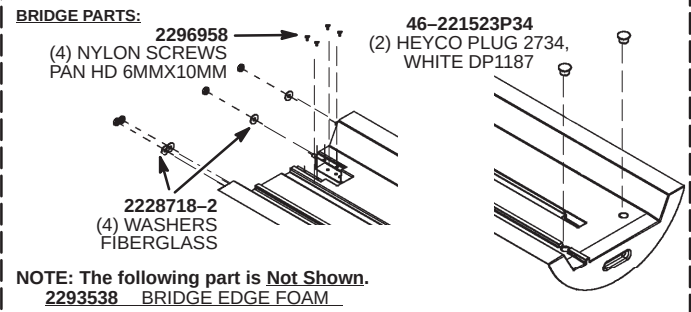LOW PROFILE CARRIAGE COVER PARTS:**BRIDGE PARTS:**


8 MAGNET AND ENCLOSURE (ASSEMBLIES)

TOUCH-UP PAINTS:

2183491-3	GRAY #1	0.6 OZ. NON-AEROSOL
46-303460P3	STEEL BLUE	0.6 OZ. NON-AEROSOL
46-303460P4	DARK STEEL BLUE	0.6 OZ. NON-AEROSOL

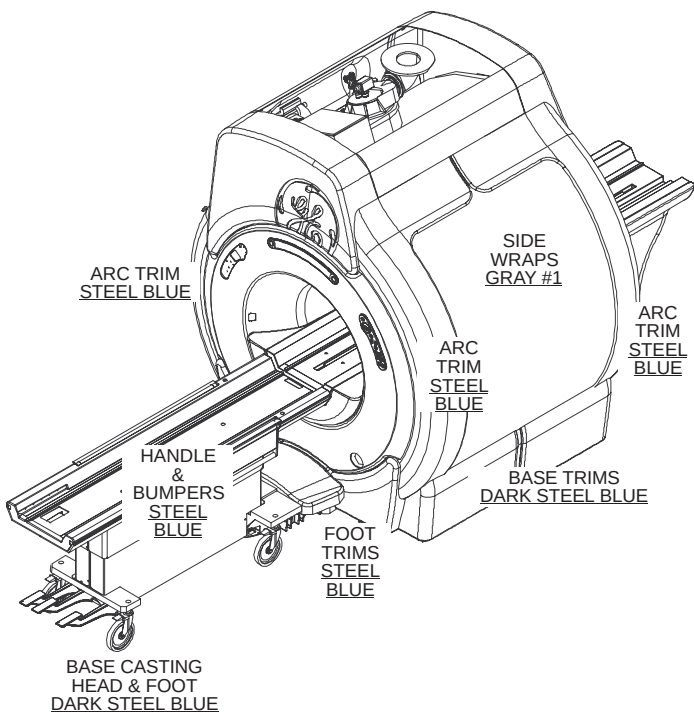





NOTE: To obtain enclosure part numbers not found in this manual, refer to Direction 2138247, *Signa Illustrated Parts List*, for enclosure covers, endbells, and other parts.

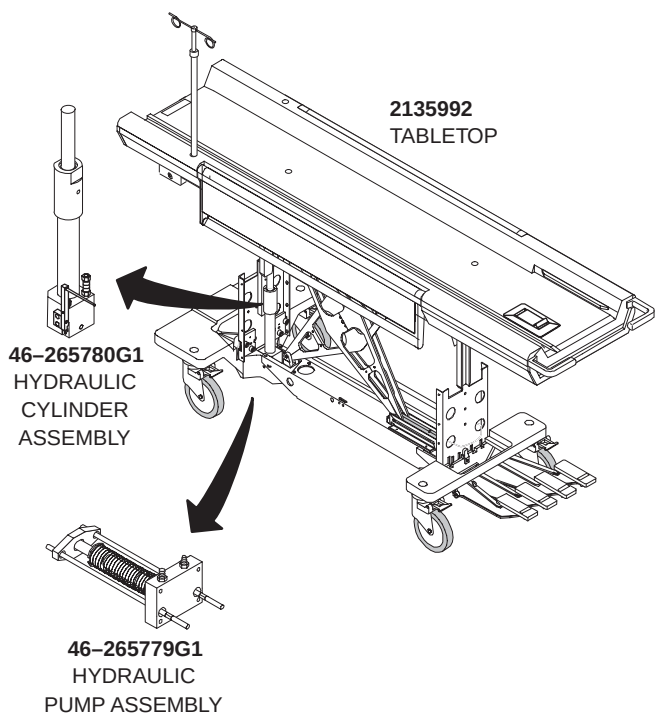


[illegible]

SYSTEM COOLING CABINETS:


REFER TO ELLIS & WATTS SERVICE MANUALS FOR FRU INFORMATION

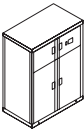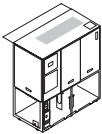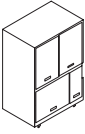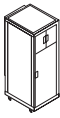OUTDOOR TSCC   REMOTE CONTROL PANEL	INDOOR WATER COOLED TSCC 	INDOOR AIR COOLED TSCC 
OUTDOOR TGWC   REMOTE CONTROL PANEL	INDOOR WATER COOLED TGWC   REMOTE CONTROL PANEL	INDOOR AIR COOLED TGWC 

REFER TO AIRSYS SERVICE MANUALS FOR FRU INFORMATION

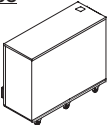
GWHX 	MRCC   REMOTE CONTROL PANEL
--	---

2343325 AIRSYS HEAT EXCHANGER (GWHX)**2341077-2** AIRSYS QUICK DISCONNECT KIT**2341077-3** AIRSYS REMOTE CONTROL PANEL**2298100-4** E & W MODEL LTL10 HEAT EXCHANGER CABINET**2298100-5** CONTROL PANEL FOR LTL10 ONLY**2297672** GLYCOL USED WITH E & W OR AIRSYS

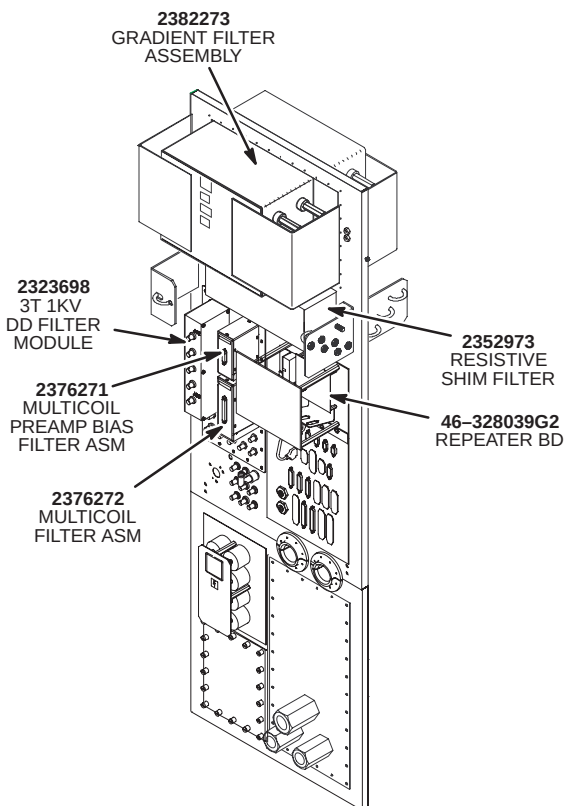


14

NOTES

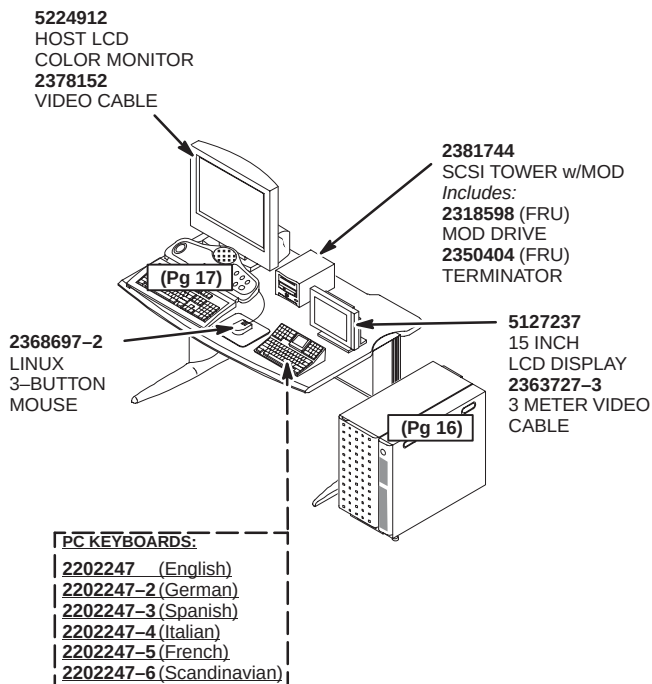


NOTE: The following parts are Not Shown.

2156158-46	S/W COLLECTOR EXCITE 11.0 G3 M4A
2389540	SCAN TOOLS MOD

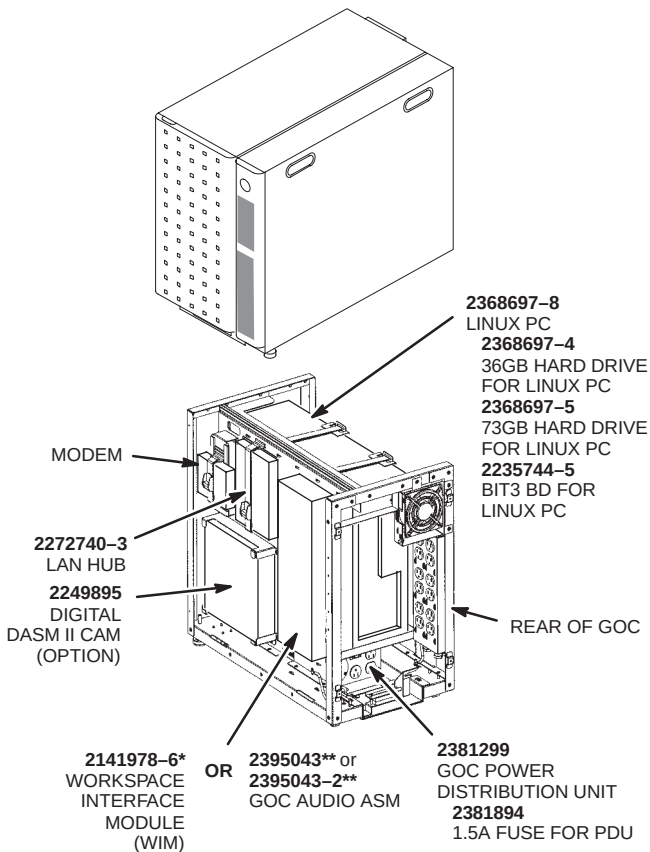GOC COMPUTER CABINET – LINUX:**NOTES:**

* If WIM is not available, order part 2141978-5 and replace cable between motherboard J1 and CIF J3.

** If replacing WIM with GOC AA, order cable 2395735.

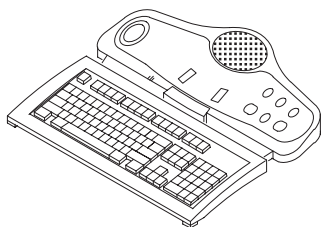LATEST VERSIONS:**KEYBOARDS ONLY:**

2275756 (English)
*2275756-2 (French)
*2275756-3 (German)
*2275756-4 (Scandinavian)
*2275756-5 (Italian)
*2275756-6 (Portuguese)
*2275756-7 (Spanish)

*** NOTE:**

These types of Keyboard/SCIM Assemblies were not available at the time this book was released, see below for earlier types.

NOTE: The following parts are Not Shown.

2262158-9	CABLE, 50 PIN SCIM KEYBOARD
2289697-2	SCIM ^W /ENGLISH OVERLAY
2289697-8	MR SCIM GLOBAL OVERLAY SET

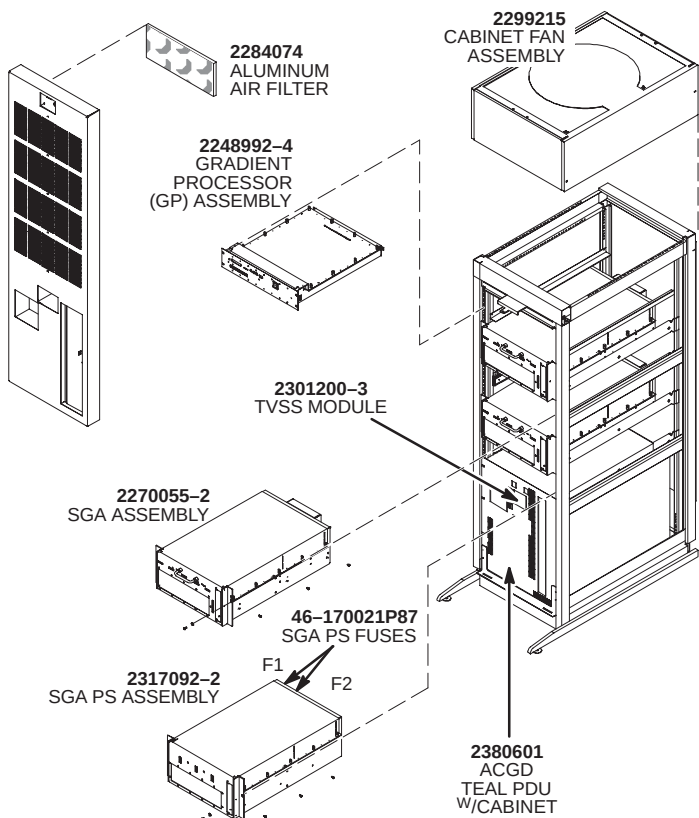
NOTE: Refer to the Directions on CD–2160623 (Signa Service Methods) for additional product information.



TABLE 1–1

LOCATION	AMPS	VOLTS	GEMS PART NUMBER	MANUFACTURER	
				BUSS	G–S
FUSE F1	1.25	480	2276547–2	SNQ–R–1–1/4	ATQR1–1/4
FUSE F2	3	120	2276547–3	FNQ–R–3	ATQR3
FUSE F3	1	12VDC	2276547–8	LP–CC–1	ATDR1
FUSE F4	0.12	480	2276547–9	SNQ–R–8	ATQR1/8
FUSE F5	125	480		LPJ125	AJT125
FUSE F6	50	480		LPJ–50	AJT50
FUSE F7	20	480	2276547–7	LPJ–20	AJT20
FUSE F8	5	480		FRS–R–5	TRS5R
FUSE F9	8	400	2276547–6	FRS–R–8	TRS8R
FUSE F10	10	120	2276547–11	FNQ–R–10	ATQR10
FUSE F11	10	120	2276547–11	FNQ–R–10	ATQR10

2276547–14	12 VDC BATTERY
2276547–16	12 VDC SUPPLY

ACGD/PDU CABINET ASSEMBLY W/SGA/SGAPS/GP


NOTE: The following Part is Not Shown.

2261124-4 ACGD GRADIENT DATA PDL DOCUMENT

2374930-2 3T ACGD PDU CONTROL BOARD

ACGD/PDU CABINET

[illegible]

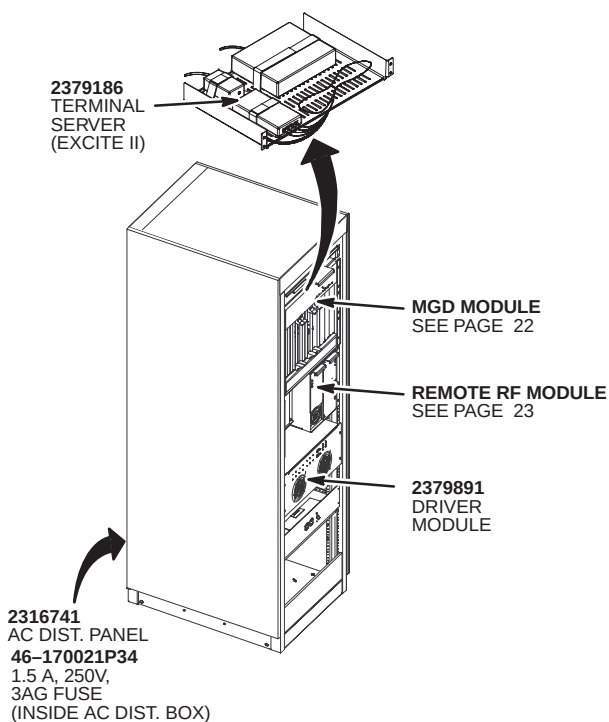


NOTE: The following System Cabinet Options are Not Shown.

2325947-12 MINERVA AP (MASTER)

2325947-13 MINERVA AP (SLAVE)

2325947-14 MINERVA PIM


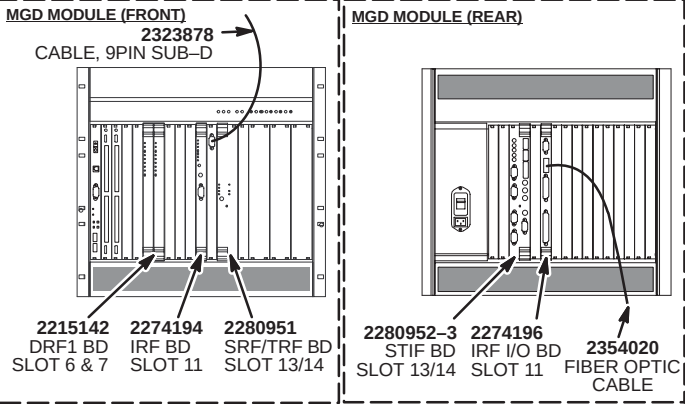


NOTE: The following MGD Module parts are Not Shown.

2294300-3	AGP BOARD
2294300-4	APS BOARD WITH IT DAUGHTER BD
2294300-5	BRIDGE BOARD
2294300-6	MGD CHASSIS W/BACKPLANE & AC ENTRY (No Boards)
2294300-7	MGD FAN
2294300-8 or	
5181495-3	POWER SUPPLIES
2294300-9	INDICATOR BOARD
2294300-19	SCP2 BOARD

NOTE: The following RRF MODULE (page 23) Parts/Cables are Not Shown.

46-306565P1	CABINET FILTER
2308729	BNC TO SMB CABLE (MUX J8 TO CAB I/F J1)
2308728	8W8 CABLE ASSEMBLY (MUX J5)
2322578	J15/23/25/26 DATA INPUT CABLE (CAB I/F)
2322580	CABLE J18, UNBLANK (CAB I/F)
2322581	CABLE J24, MAGNET MONITOR (CAB I/F)
2322582	CABLE J30/31, SERIAL EXP (CAB I/F)
2322583	CABLE J62 (15 PIN SUB-D) (CAB I/F)
2322586	CABLE J61 (37 PIN SUB-D) (CAB I/F)
2315093	ETHERNET CABLE
2315093-2	CABLE J63 ETHERNET (CAB I/F)
2308729-2	BNC TO SMB CABLE (MUXJ13 TO CAB I/F J2)
2308729-3	BNC TO SMB CABLE (MUXJ14 TO CAB I/F J3)

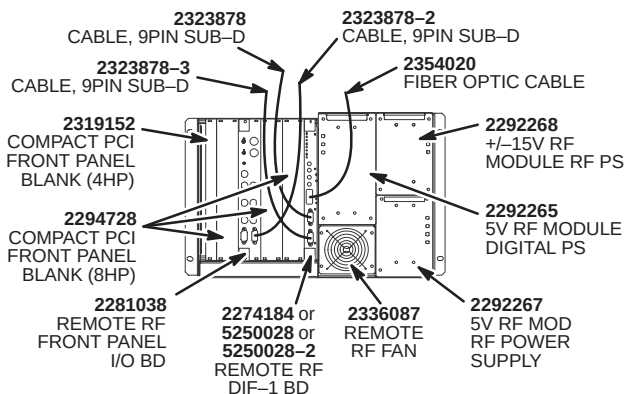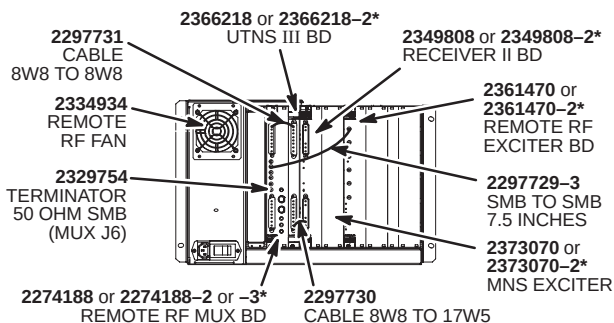
2305884 REMOTE RF MOD ASM (CHASSIS W/MIDPLANE) (FRONT)**2305884 REMOTE RF MOD ASM (CHASSIS W/MIDPLANE) (REAR)**

NOTE: Additional REMOTE RF MOD ASM Parts/Cables are located on page 22.

* A part that does not have a "-2" or "-3" can be replaced with either the same part number or the same number ending with a "-2" or "-3". If replacing the "-2" or "-3" version, then only the part number with a "-2" or "-3" can be used.

24

NOTES

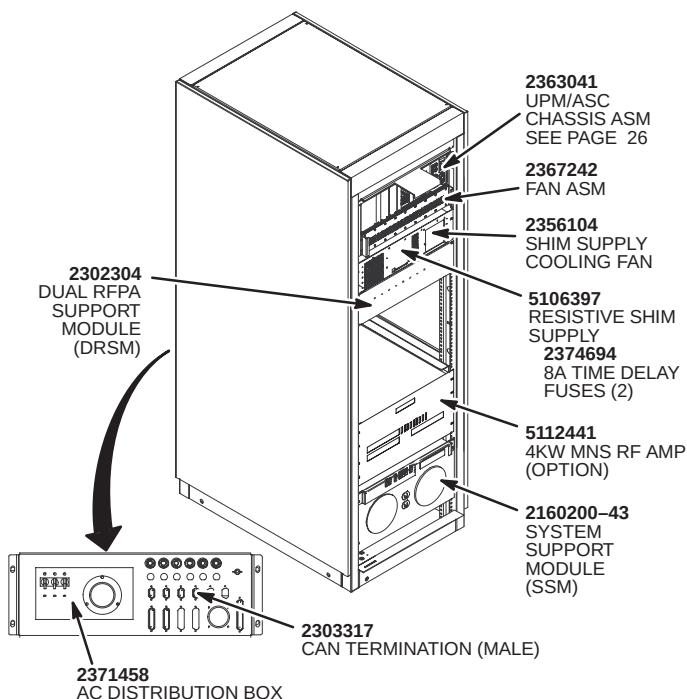

NOTE: The following RF ACCESSORY CABINET parts Not Shown:

46-306565P1 10" FRONT COVER FILTER

2365797 CABLE - I/O J127 TO ASC SLOT 11/14 J1/2/4/5

2368552 CABLE - ASC SLOT 14-J3 TO I/O J100

2369552-8 CABLE - ASC SLOT 13-J4 TO I/O J102

2369552-9 CABLE - ASC SLOT 13-J5 TO I/O J103

2369552-10 CABLE - ASC SLOT 10-J4 TO I/O J104

2369552-11 CABLE - ASC SLOT 10-J5 TO I/O J105

2379507 DRSM JUMPER PLUG

2383066 DRSM JUMPER CABLE (ASC-DRSM-I/O)

FOR 2363624 RF ACCESSORY CABINET:

2363041 or 2363041-2*
ASC CHASSIS
ASSEMBLY


2359516 or 2359516-2*
UPM RF DETECTOR
BOARD

2363040 or 5120955
ASC CHASSIS
POWER SUPPLY

2329754
50 OHM SMB
TERMINATOR (2)

2359518
UPM PROCESSOR
BOARD

2369923
ASC INPUT PANEL

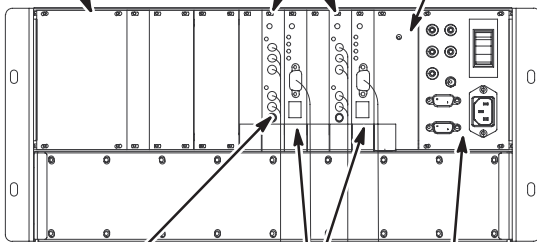**FOR 2363624-2 RF ACCESSORY CABINET WITH MNS:**

2363041 or 2363041-2*
ASC CHASSIS
ASSEMBLY


2359516 or 2359516-2*
UPM RF DETECTOR
BOARD

2363040 or 5120995
ASC CHASSIS
POWER SUPPLY


2329754
50 OHM SMB
TERMINATOR (2)

2359518
UPM PROCESSOR
BOARD

2369923
ASC INPUT PANEL

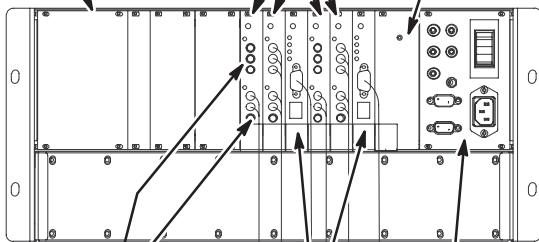

* When replacing RF Detector boards or Chassis, the "no dash" part can only be used when replacing a "no dash" part. If the part is a "-2", then the replacement part must be a "-2".

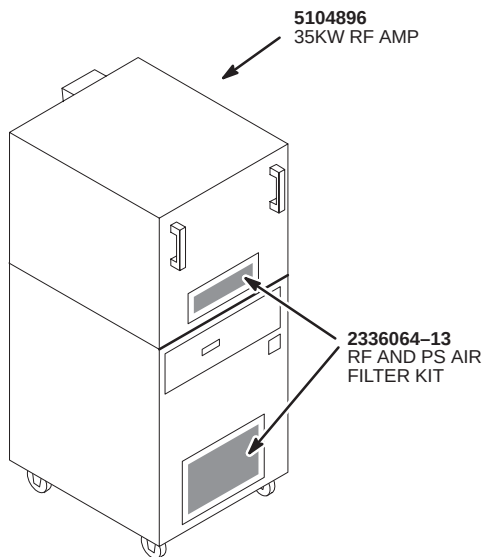

NOTE: The following NB RF AMP CABINET parts are Not Shown:

2336064-8 35KW RF AMP FUSE KIT

2336064-15 35KW RF AMP TEMPERATURE SENSOR

NOTE: The following RF TRANSMIT FLEX CABLES are Not Shown:


2384850 RUN 1156 – MR8-J4 TO CONNECT TO RUN 1125

2384850-2 RUN 1160 – MG3-A3-J1 TO CONNECT TO RUN 1127

2384851 RUN 1157 – MR8-J3 TO CONNECT TO RUN 1126


2384851-2 RUN 1158 – MR1-CPC-J2 TO CONNECT TO RUN 1137

2384851-3 RUN 1159 – MR9-J4 TO CONNECT TO RUN 1146

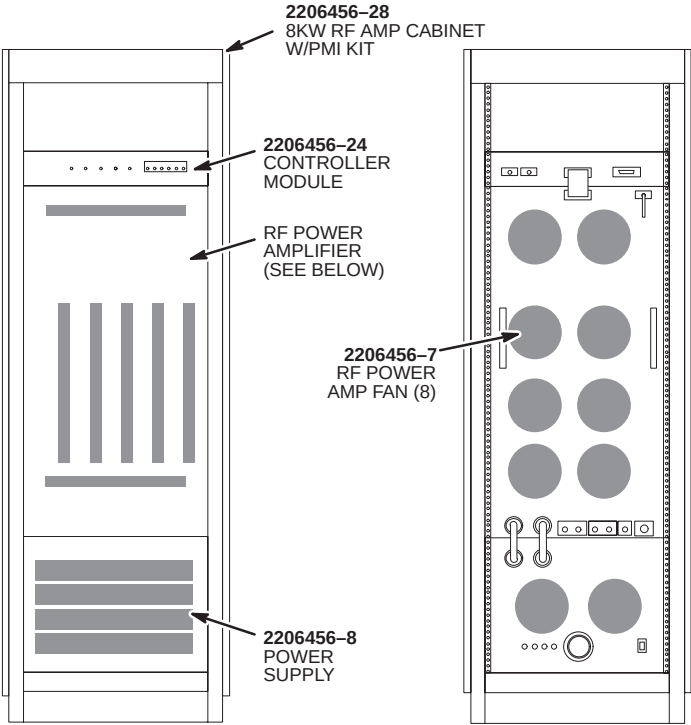


NOTE: The following RF POWER AMP internal parts are Not Shown:

2206456-2	RF POWER AMP PREAMP/SPLITTER
2206456-3	RF AMPLIFIER (5)
2206456-4	RF POWER AMP COMBINER/COUPLER
2206456-5	RF POWER AMP COUPLER (4)
2206456-6	RF POWER AMP BIAS TEE
2206456-13	CABLE – SPLITTER TO AMPLIFIER (5)
2206456-20	COUPLER-BIAS TEE
2206456-21	POWER BIAS TEE
2206456-22	PMI INTERFACE KIT
46-170021P29	5 AMP FUSES

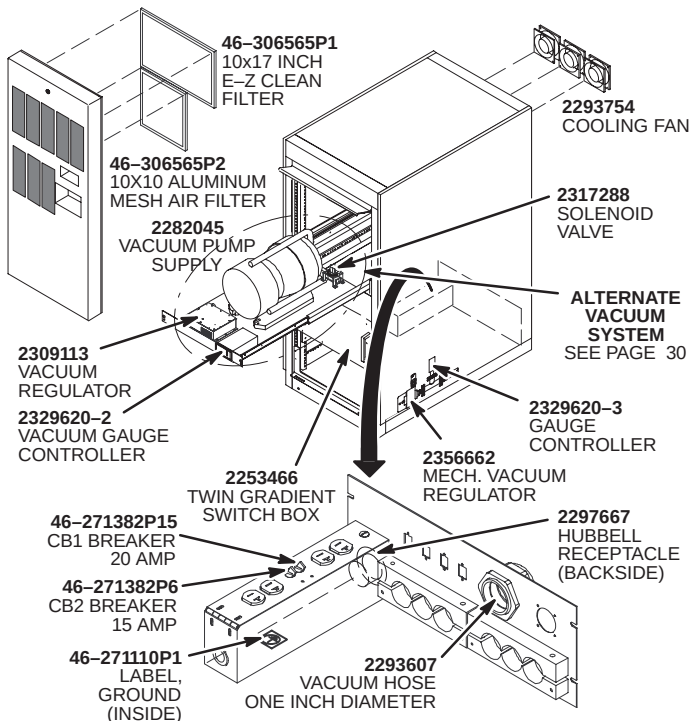

NOTE: The following parts are Not Shown.


2282045-2	VACUUM PUMP TIP SEAL KIT
2293608	QF25 CENTERING RING & O-RING
2293608-2	QF16 CENTERING RING & O-RING
2293609	QF25 VACUUM ELBOW
2293616	QF25 ADAPTER TO ONE INCH HOSE
2293617	VACUUM REDUCING 'T'
2293618	QF16 CLAMP
2293618-2	QF25 CHAIN CLAMP
2294010-3	HOSE CLAMP, SST ONE INCH
2297760	VACUUM GREASE, MSDS# 8366442
2329620-4	VACUUM GAUGE CONTROLLER CABLE
2284301	VACUUM GAUGE TO TAC J2 CABLE

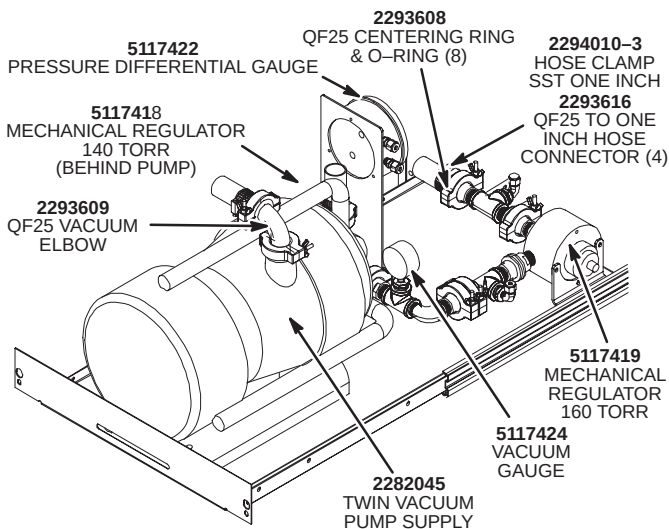

NOTE: The following parts are Not Shown.

2282045-2	VACUUM PUMP TIP SEAL KIT
5118818	CLAMP NW 20-25
2297760	VACUUM GREASE, MSDS# 8366442


5127147
SERVICE KIT FOR TAC VACUUM
SYSTEM UPGRADE
(Used to replace old style vacuum
system, with solenoid)

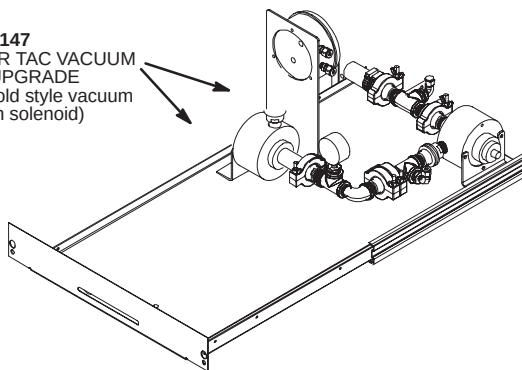

REPLACEMENT CABLES

The table below provides information for generic replacement cables. If any of the cables are ordered, make sure to order write-on labels to be used for identifying cable part number and connection designations at the end of each cable.

Description	Length	
	65 Feet	100 Feet
9 Pin Sub-D Male/Female	46-244714G30	46-244714G37
9 Pin Sub-D Female/Female	46-244713G30	46-244713G37
15 Pin Sub-D Male/Female	46-244712G30	46-244712G37
15 Pin Sub-D Female/Female	46-244711G30	46-244711G37
25 Pin Sub-D Male/Female	46-244710G30	46-244710G37
25 Pin Sub-D Female/Female	46-244709G30	46-244709G37
37 Pin Sub-D Male/Female	46-258877G30	46-258877G37
	50 Feet	65 Feet
BNC RG 223	46-251710G1	46-251710G18
Cable Write-on Labels (60 labels per booklet)	46-208781P3	

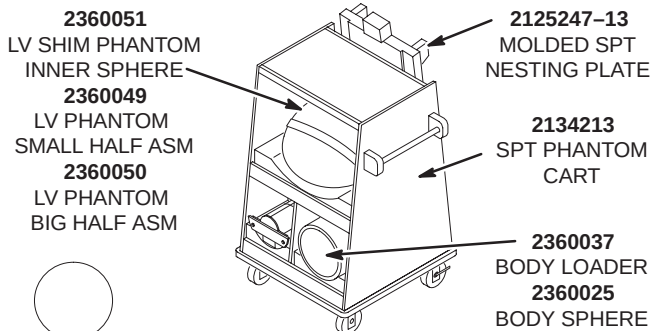


2359877
HEAD
SPHERE



2360034
100mm
SPHERE


5121115
LOCATION
BRACKET
w/SCREWS

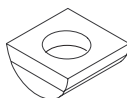


2291292
EPI HEAD
POSITIONER




2170481
EPI FOAM
POSITIONER

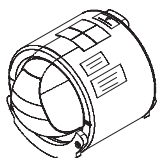

2360031
HEAD LOADER

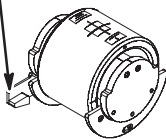


2371512
DQA PHANTOM
WITH LOADER

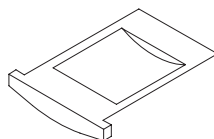


5110241
HEAD LOADER
POSITIONER



2152220
MRS PHANTOM

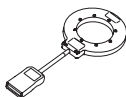

2386800
MNS SERVICE COIL

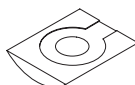

2389665
MNS PHANTOM
POSITIONER

[illegible]

ACCESSORIES:

The following E-Catalog accessories are available for purchase with M-Catalog Options.

COLOR POSTSCRIPT PRINTER FOR A/W E7014LA

2112346	POSTSCRIPT COLOR NETWORK PRINTER
----------------	----------------------------------

PATIENT ECG LEADWIRE SET E8811EC

2219452	HIGH IMPEDENCE ECG LEAD KIT
46-306867G1	PATIENT ECG CABLE ASM (E8811EC)

COILS:**3T 8 CHANNEL TORSO ARRAY COIL**

2381683-2	COIL WITH CABLE
2381683-3	PHANTOM
2382683-4	PHANTOM POSITIONER

3T 8 CHANNEL BRAIN COIL

2380637-2	BRAIN COIL/MIRROR/BASEPLATE ASM
2380637-3	PHANTOM POSITIONER
2380637-5	CABLE ASM

3T GP FLEX COIL

2382817	GP FLEX COIL
2389979	GP FLEX COIL CABLE KIT

3T QUAD EXTREMITY COIL

2380635-2	KNEE/FOOT ASM W/CABLE
2380635-5	KNEE/FOOT CABLE ASM
2357819-10	EXTREMITY REPAIR KIT
2357819-11	PHANTOM POSITIONER
2357819-12	PHANTOM

3T T/R EXTREMITY COIL

2357819-2	3T KNEE/FOOT ASM
2357819-9	3T CABLE ASM
2357819-10	EXTREMITY REPAIR KIT
2357819-11	PHANTOM POSITIONER
2357819-12	PHANTOM

3T SPLIT TOP HEAD COIL

2376114	SPLIT HEAD COIL ASM
2376104	QD BOX ASM
2145391	MIRROR
2355094	REPLACEMENT SPRING CONTACTS KIT

3T NV ARRAY COIL

2414383	E3 NV ARRAY W/O CABLE
2414385	E3 NV SYSTEM COIL
2414386	ANTERIOR CHEST WITH CHEST CABLE
2414389	NV PHANTOM SET
2414391	E2 NV SYSTEM CABLE

OPTION KEYS:

2383494	Fiesta C for EXCITE
2390037	Asset for EXCITE
2390040	ConnectPro for EXCITE
2390075	3D FRFSE for EXCITE
2390078	EchoPlus for EXCITE
2390079	Diffusion Tensor EXCITE
2390084	Probe 2000 for EXCITE
2390085	Probe 200 Upgrade for EXCITE
2390086	Sage 7 for EXCITE
2390091	2D Delayed Enhancement EXCITE
2390095	iDrive ProPlus for EXCITE
2390097	SmartPrep 2000 for EXCITE
2390098	Flouro Trigger MRA for EXCITE
2390103	Performed Proc Step for EXCITE
5110228	Single Voxel Probe–P for EXCITE
5110229	Single Voxel Probe–S and Probe–P for EXCITE
5110230	Probe 2D – CSI for EXCITE
5110231	Probe 3D – CSI for EXCITE
5110232	MNS, SAGE, Probe–P for EXCITE
5110319	Research Proton Spectro for EXCITE

BrainWave OPTIONS:

2363727	15" LCD Display with 5 foot VGA Cable for Computer connection
2394034-5	(Cedrus) Response Pad Asm. Left & right hand units with two buttons for each hand. Includes fiber optic cable (20 ft) and optical converter box
2394034-6	(Cedrus) "Lumina" response controller box
2394034-7	(Cedrus) AC Power module for Lumina Controller Box
2394034-9	FRU usage only Kit containing all cables [Cedrus-4, -5, -6, -7, 2xxxxxx-6, 46-265505P1 (qty 2)]
2394034-10	Penetration Sub-Panel Asm
2394034-16	PC compatible Keyboard (narrow)
2394034-17	PC compatible Mouse
5110767	BrainWave Application Software CD-ROM
5110736	EXCITE II Scanner License Key MOD for the BrainWave H/W Option
5110946	PC and Platform Software Asm
5110968	Patch CD for BrainWave Stimulus PC