

Intera 0.5 T		Intera 1.0 T		Intera 1.5 T	
Standard	Omni	Power	Omni	Power	Master

Section 13:

Covers

Contents


1.	MAGNET COVERS	3
1.1.	Prepare magnet covers	4
1.1.1.	PICU and Physiology unit left or right	4
1.1.2.	Cut out the marked area of side cover.....	4
1.1.3.	Install the accessory hooks	5
1.1.4.	front cone cover adaption.....	6
1.2.	Covers installation	7
1.2.1.	Install front and rear lower cover.....	7
1.2.2.	Install cone cover front	8
1.2.3.	Install the top covers	9
1.2.4.	Install the bridge section	10
1.2.5.	Install the side covers.....	11
1.2.6.	Install the PICU cover front	12
1.2.7.	Install the PICU cover rear	12
2.	CHECK AIR COOLING SYSTEM	13
3.	MOUNTING THE PATIENT SUPPORT INTERA COVERS.....	14
4.	MOUNTING THE SYSTEM FILTER BOX COVERS.....	15

Intera 0.5 T		Intera 1.0 T		Intera 1.5 T	
Standard	Omni	Power	Omni	Power	Master

Intera 0.5 T		Intera 1.0 T		Intera 1.5 T	
Standard		Omni	Power	Omni	Power Master

1. MAGNET COVERS

Figure 1 – Magnet covers overview

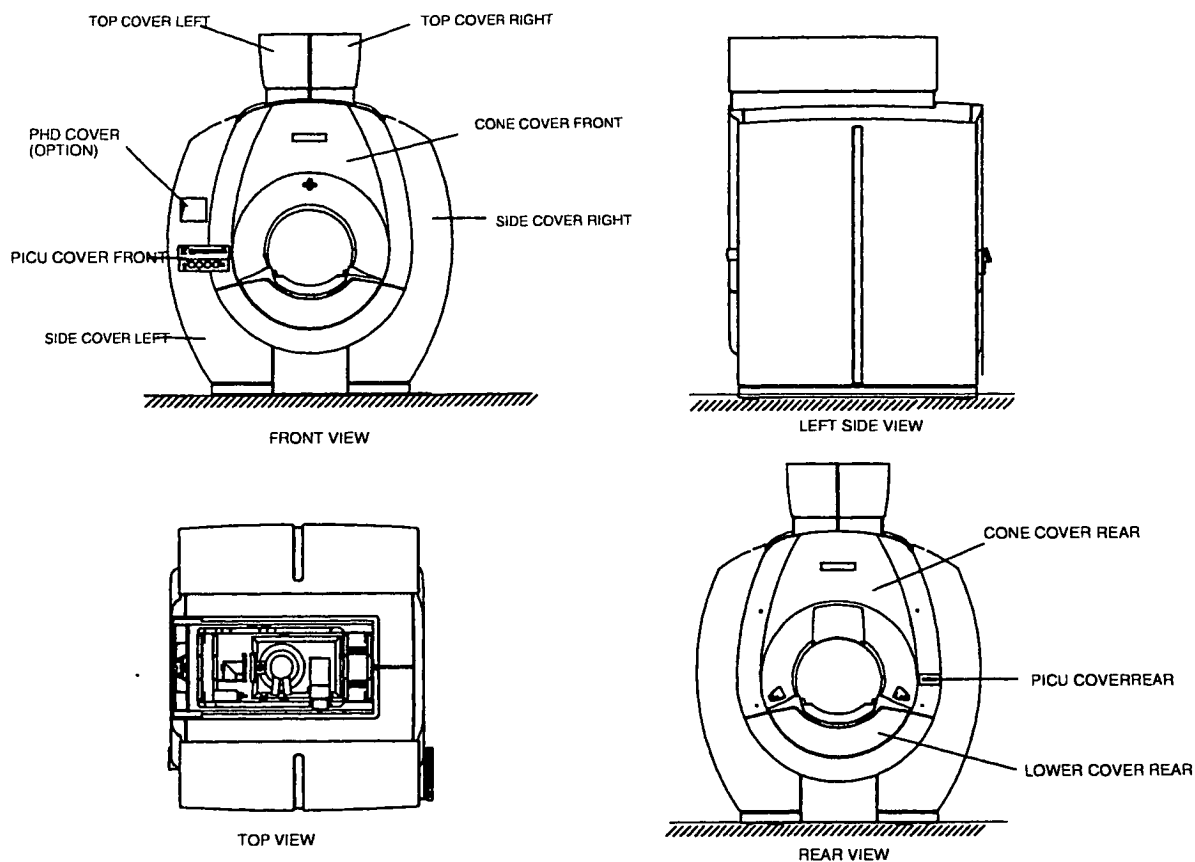

Intera 0.5 T		Intera 1.0 T		Intera 1.5 T	
Standard	Omni	Power	Omni	Power	Master

1.1. PREPARE MAGNET COVERS

1.1.1. PICU AND PHYSIOLOGY UNIT LEFT OR RIGHT

The position of the PICU (on the magnet) and of the Physiology unit (underneath the patient support) must be discussed with the customer. The next two points must be taken into account:


- The customer might want the PICU at the same side of the system where the entrance door is.
- The Physiology unit must be on the same side as the PICU.
- The Trolley (option) can only be used at the opposite side of the PICU and Physiology unit.
- Position of the coil cabinet.

Note that the head coil must be modified in case the PICU is mounted on the right side.

1.1.2. CUT OUT THE MARKED AREA OF SIDE COVER

Cut out the marked area in the side cover right or left. This can be done using an electric jigsaw.

Figure 2 – Prepare side covers

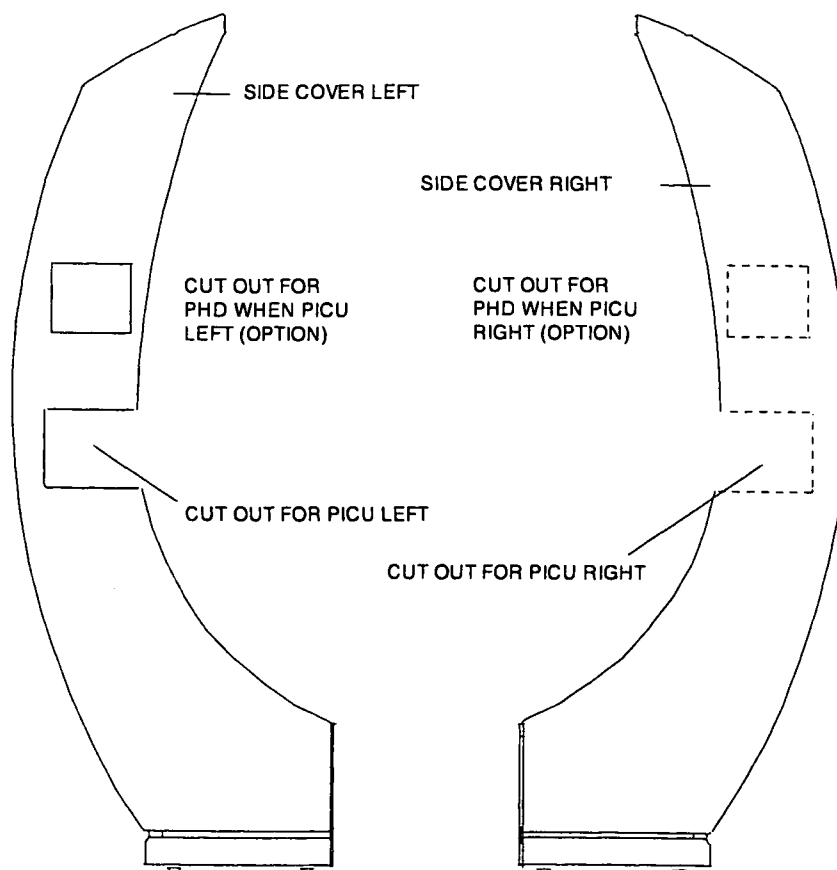

FH81.DRW


If a PHD (option) must be installed, also cut out the marked area above the PICU.

Intera 0.5 T		Intera 1.0 T		Intera 1.5 T	
Standard	Omni	Power	Omni	Power	Master

1.1.3. INSTALL THE ACCESSORY HOOKS

Install the accessory hooks to the side cover below the location of the PICU area.
Two markings are present at the inside of the side cover.
Drill holes through the cover, and mount the accessory hooks.

Figure 3 - Accessory hooks

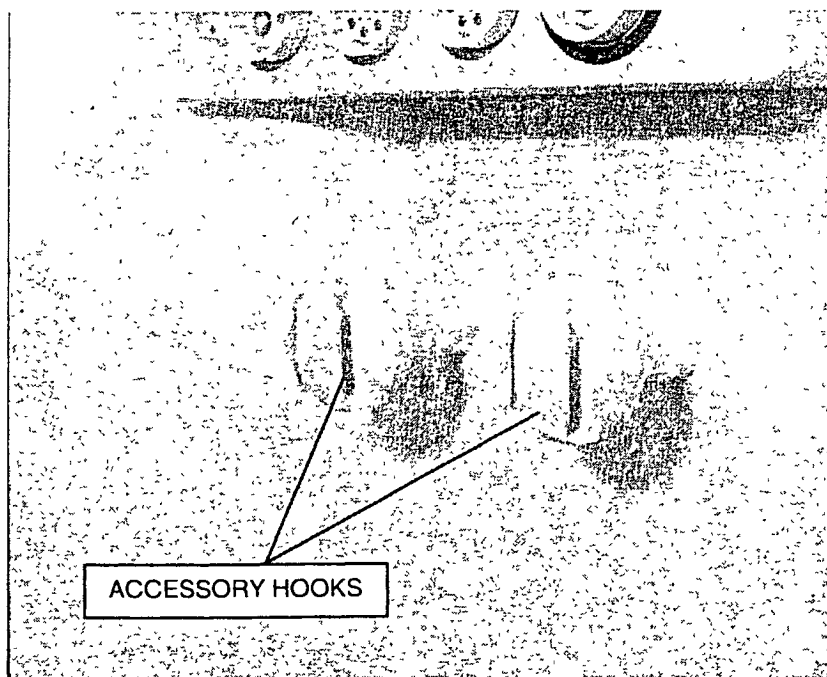


Intera 0.5 T		Intera 1.0 T		Intera 1.5 T	
Standard	Omni	Power	Omni	Power	Master

1.1.4. FRONT CONE COVER ADAPTION

Cut out the marked area in the cone cover front. This can be done using an electric jigsaw.

Figure 4 – Cut out in cone cover front

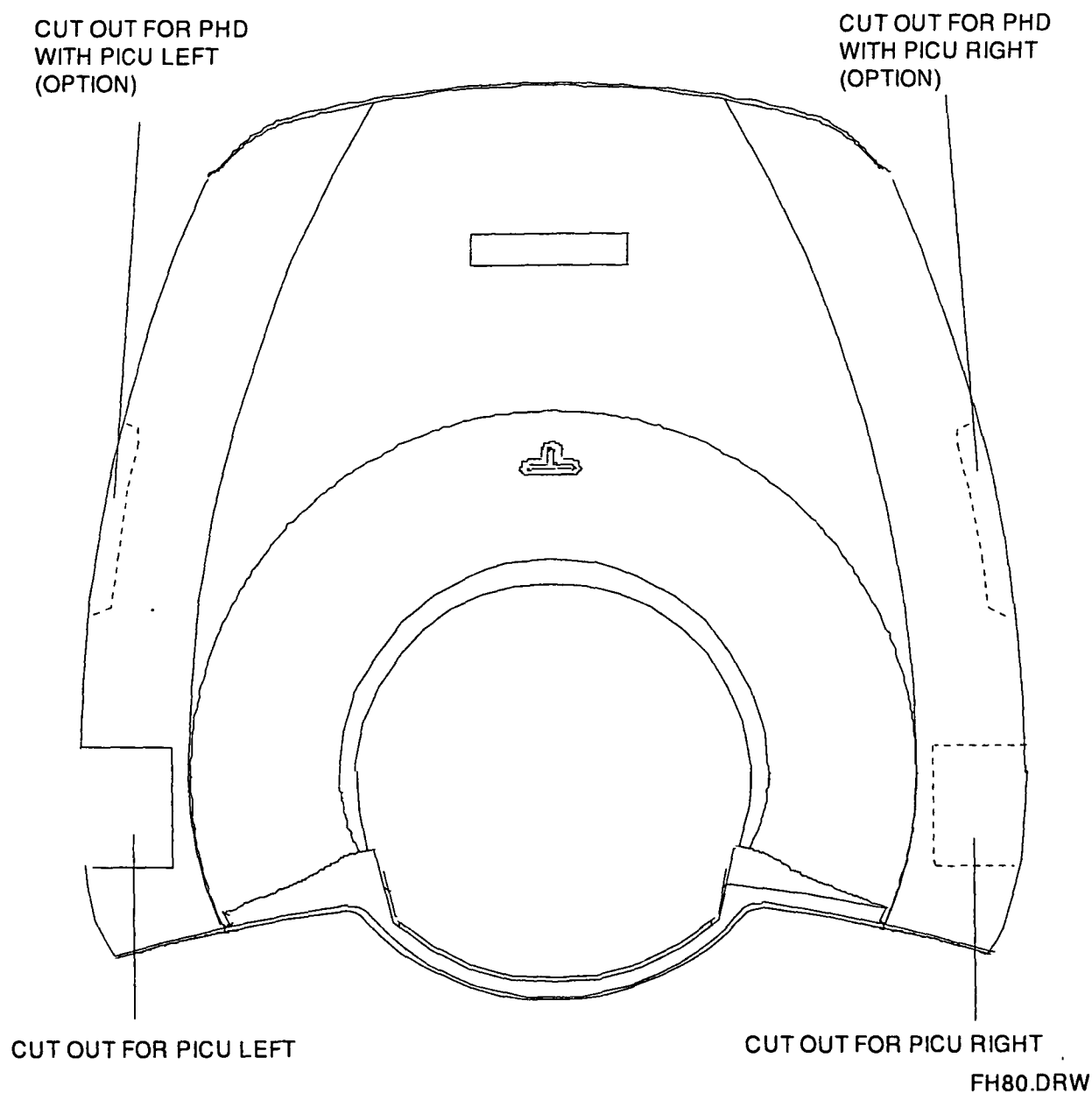

Intera 0.5 T		Intera 1.0 T		Intera 1.5 T	
Standard	Omni	Power	Omni	Power	Master

1.2. COVERS INSTALLATION

The magnet covers must be installed in the following sequence:

1. front and rear bottom covers
2. front and rear cone cover
3. top cover left and right
4. side covers
5. PICU cover front
6. PICU cover rear
7. bridge section

1.2.1. INSTALL FRONT AND REAR LOWER COVER

Figure 5 – Front lower cover


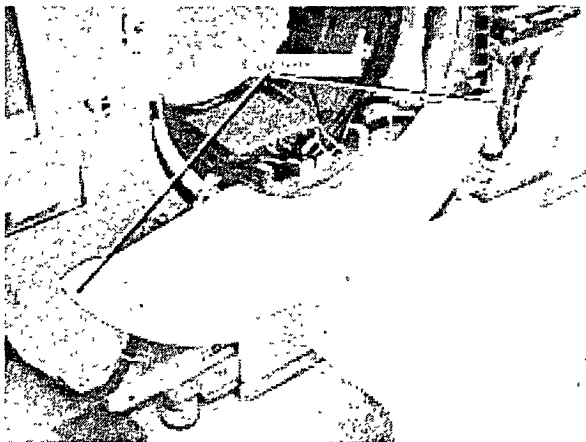


Figure 6 – Rear bottom cover

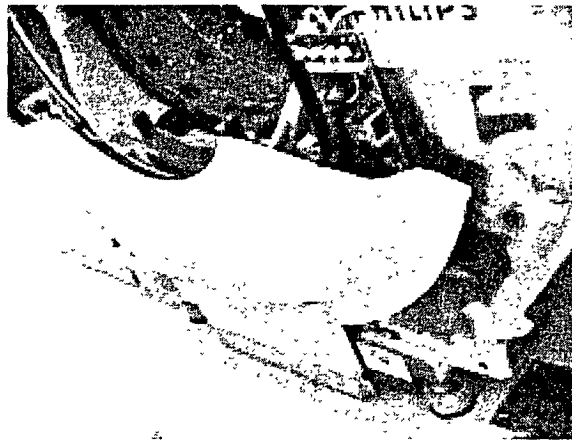

Intera 0.5 T		Intera 1.0 T		Intera 1.5 T	
Standard	Omni	Power	Omni	Power	Master

1.2.2. INSTALL CONE COVER FRONT

Figure 7 – Mount front cone cover


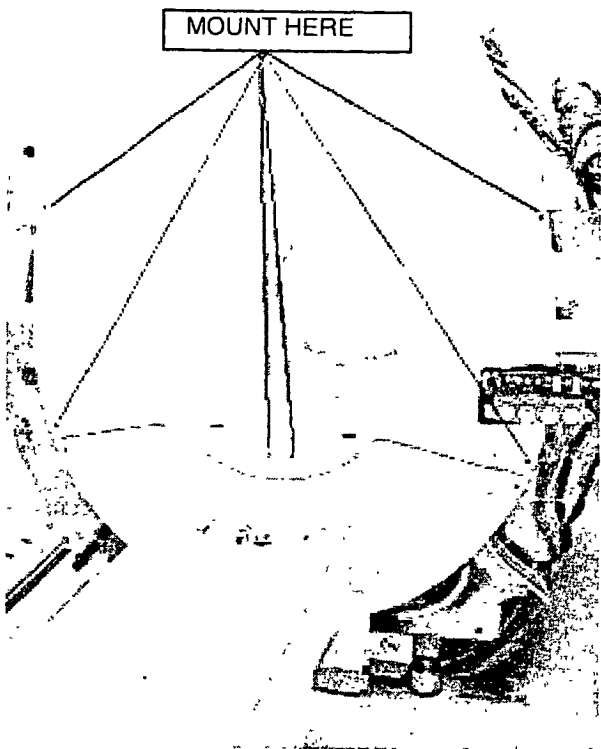


Figure 8 – Mount rear cone cover

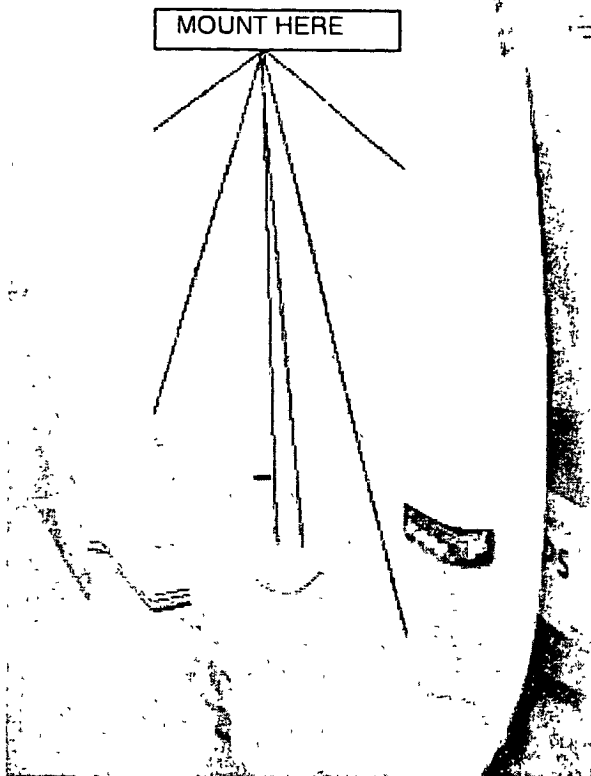

Install the rubber sealing between the QBC and the cone cover.

Intera 0.5 T		Intera 1.0 T		Intera 1.5 T	
Standard	Omni	Power	Omni	Power	Master

1.2.3. INSTALL THE TOP COVERS

1. First install the top-cover right
2. Adjust the height of the cover so that the cover is in line with the front and rear cone covers (See Figure 9, Figure 10 and Figure 11).
3. Then mount the left cone cover and adjust the height.
4. Mount the covers together as shown in Figure 12.

Figure 9 Adjust height at front bracket

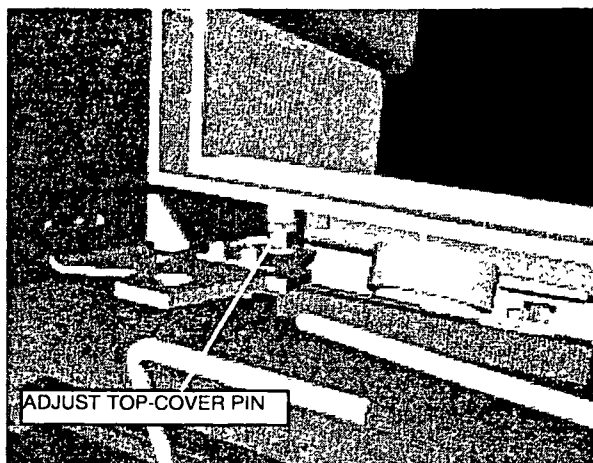


Figure 10 – Adjust height at rear bracket

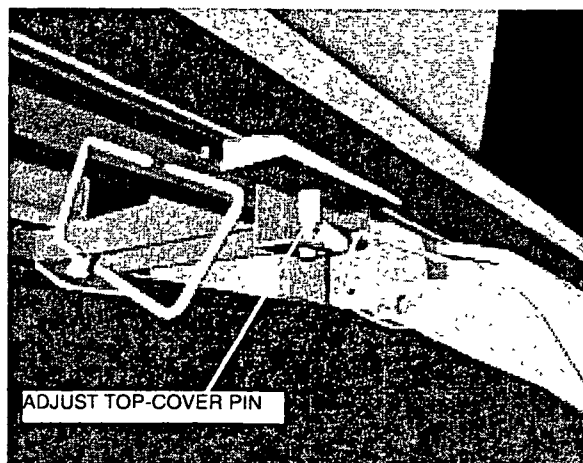


Figure 11 – Adjust height at side mounting points


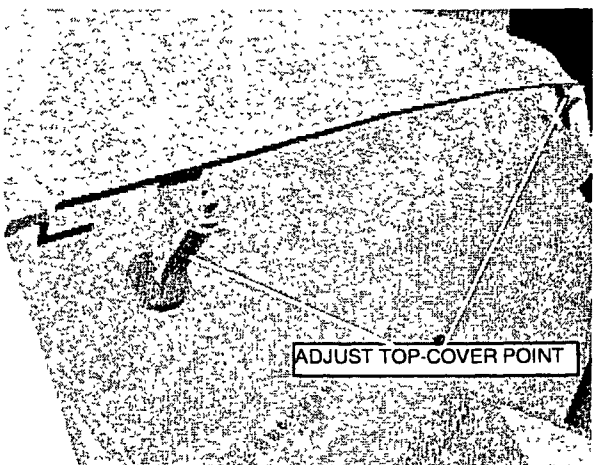


Figure 12 – Mount left and right cover together

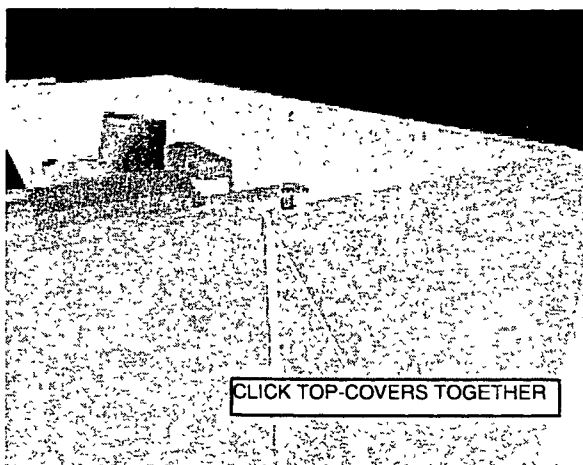

Intera 0.5 T		Intera 1.0 T		Intera 1.5 T	
Standard	Omni	Power	Omni	Power	Master

1.2.4. INSTALL THE BRIDGE SECTION

The finger protection plate is integrated into the bridge section. For cleaning purposes the finger protection plate can be removed and installed (Figure 13) by the customer.

At the bottom of the bridge section, four connection pins are present.

Position the bridge section inside the magnet bore (Figure 14), so that the pins are positioned above the four holes in the cone covers. Lower the bridge section and push it downwards.

Figure 13 – Switch plate attached to bridge section

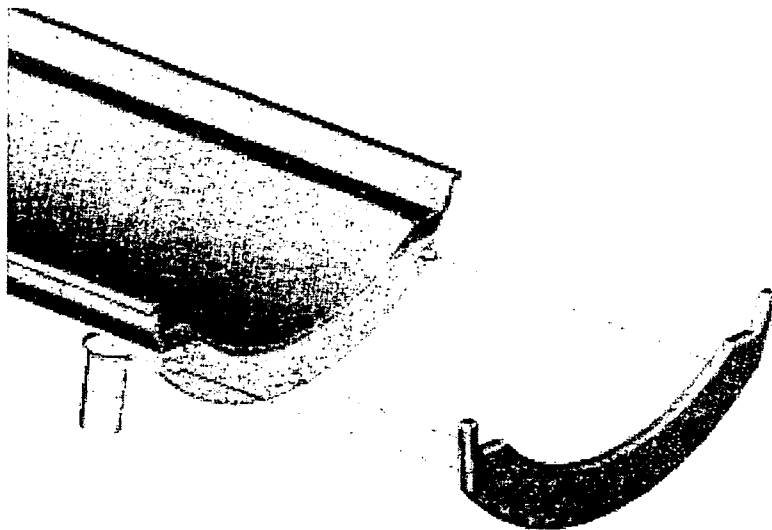


Figure 14 – Install bridge section



Intera 0.5 T		Intera 1.0 T		Intera 1.5 T	
Standard	Omni	Power	Omni	Power	Master

1.2.5. INSTALL THE SIDE COVERS

Adjust the wheels of the side cover so that they fit nicely against the cone cover and lower cover.

Figure 15 – Mount side covers

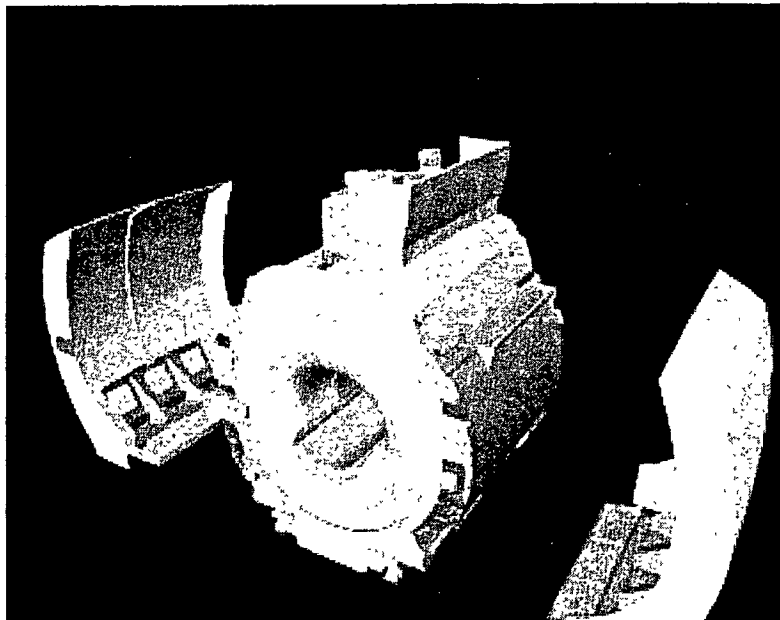


Figure 16 – Left side cover


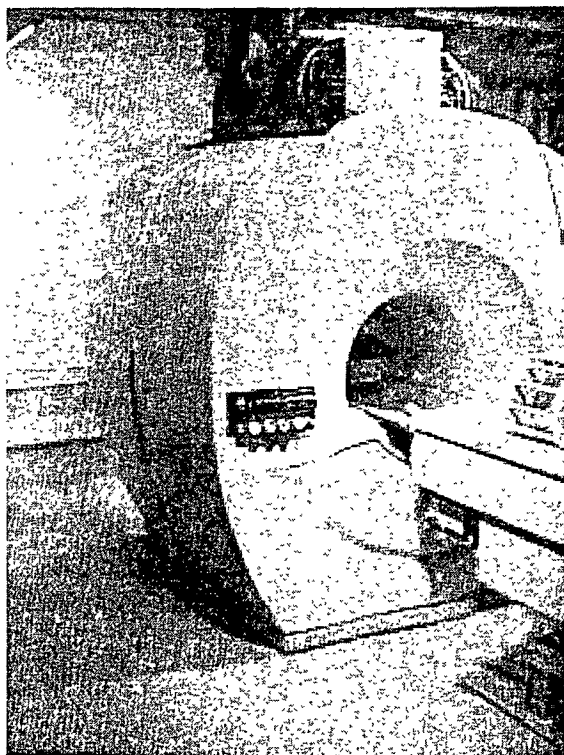


Figure 17 – Right side cover

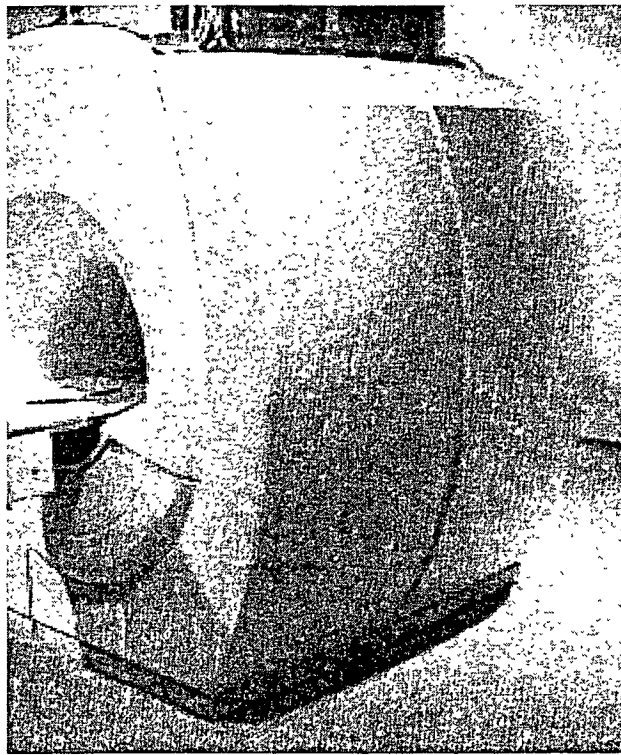

Intera 0.5 T		Intera 1.0 T		Intera 1.5 T	
Standard	Omni	Power	Omni	Power	Master

1.2.6. INSTALL THE PICU COVER FRONT

Install the PICU cover front. Install the depth of the PICU such that there is a 1-mm gap between the PICU cover and the cone cover and also between PICU cover and side cover.

Figure 18 – Install PICU cover


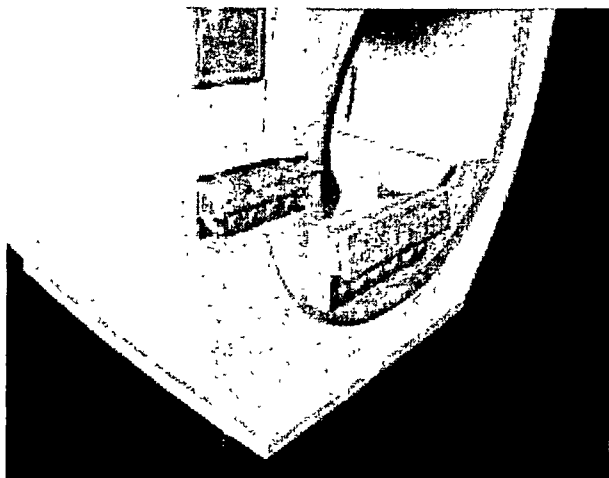


Figure 19 – Install PICU cover

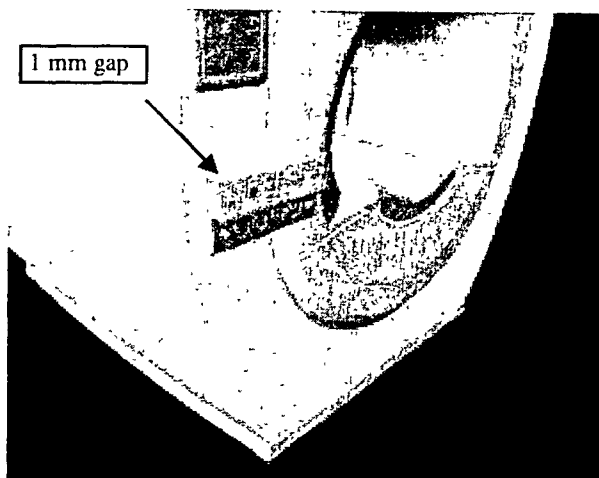


1.2.7. INSTALL THE PICU COVER REAR

Install the PICU cover rear. Install the depth of the PICU such that there is a 1-mm gap between the PICU cover and the cone.

Figure 20 – Install the rear PICU cover



Figure 21 - Install the rear PICU cover

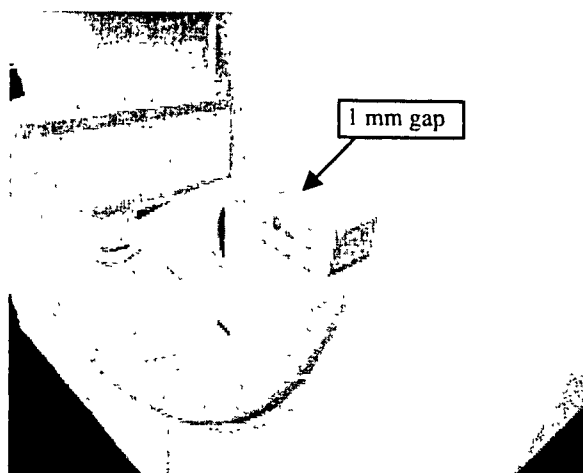

Intera 0.5 T		Intera 1.0 T		Intera 1.5 T	
Standard	Omni	Power	Omni	Power	Master

2. CHECK AIR COOLING SYSTEM

A pressure sensor is connected to one of the air ducts. This sensor activates if the pressure drops below 60 Pascal. Note that the sensor will not detect if the air volume is sufficient. The air volume can be checked by measuring the pressure at test points in the cooling system.

There are two test points where the air volume can be measured:

1. A test point is located behind the rear cone cover (See Figure 22) Measure the pressure at this test point using a manometer. The next table shows the relation between pressure [Pa] and air volume [m³/h].


Standard and Omni		Power and Master	
Air volume [m ³ /h] min-max	Pressure [Pa] Min-max	Air volume [m ³ /h] Min-max	Pressure (Pa) Min-max
800 – 1000	150 - 220	500 m ³ /h – 600 m ³ /h	95 - 125

2. Measure the pressure at the rear of the magnet ((See Figure 22 nipple on rear air duct) using a manometer.

Standard and Omni		Power and Master	
Air volume [m ³ /h] Min-max	Pressure (Pa) Min-max	Air volume [m ³ /h] Min-max	Pressure (Pa) Min-max
800 – 1000	65 - 95	500 – 600	TBD

This manometer must be arranged locally (see tool list in paragraph "Service tools).

Figure 22 – measurement point

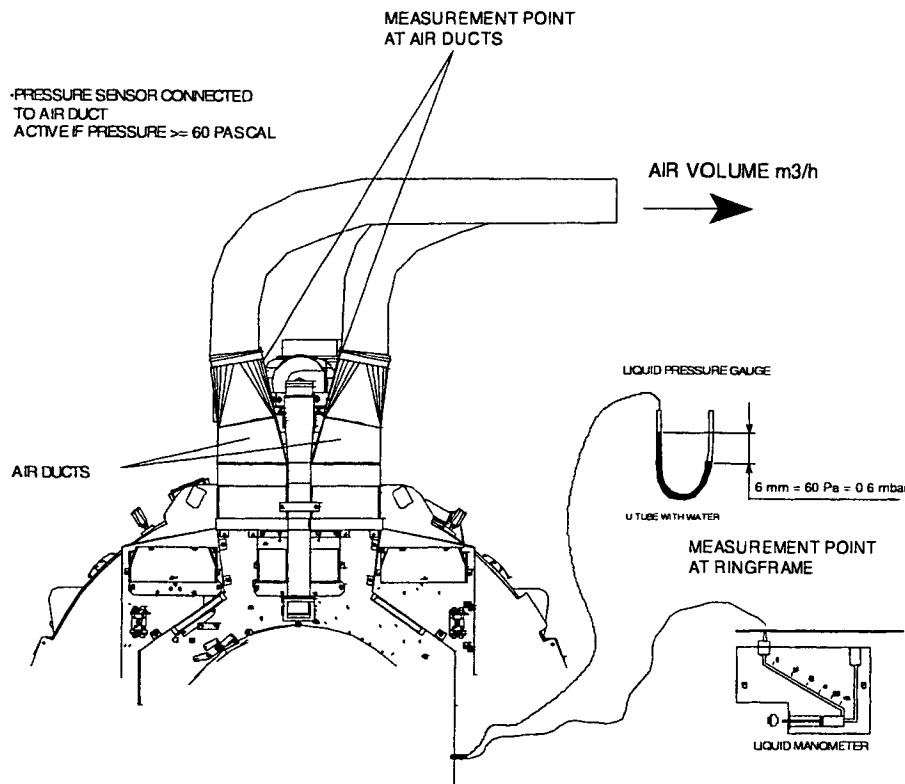


Intera 0.5 T		Intera 1.0 T		Intera 1.5 T	
Standard	Omni	Power	Omni	Power	Master

3. MOUNTING THE PATIENT SUPPORT INTERA COVERS

On the patient support Intera, most covers are factory mounted. Only the bottom covers must be mounted on site. Install the covers as follows (see also drawing Z3-5 9):

1. Lift up the middle covers.
2. Mount the left and right bottom cover.
3. Mount the earth wires to all covers (also rear and front bottom covers).
4. Mount the front and rear bottom covers.
5. Lower the middle covers

Figure 23 – Mount bottom covers

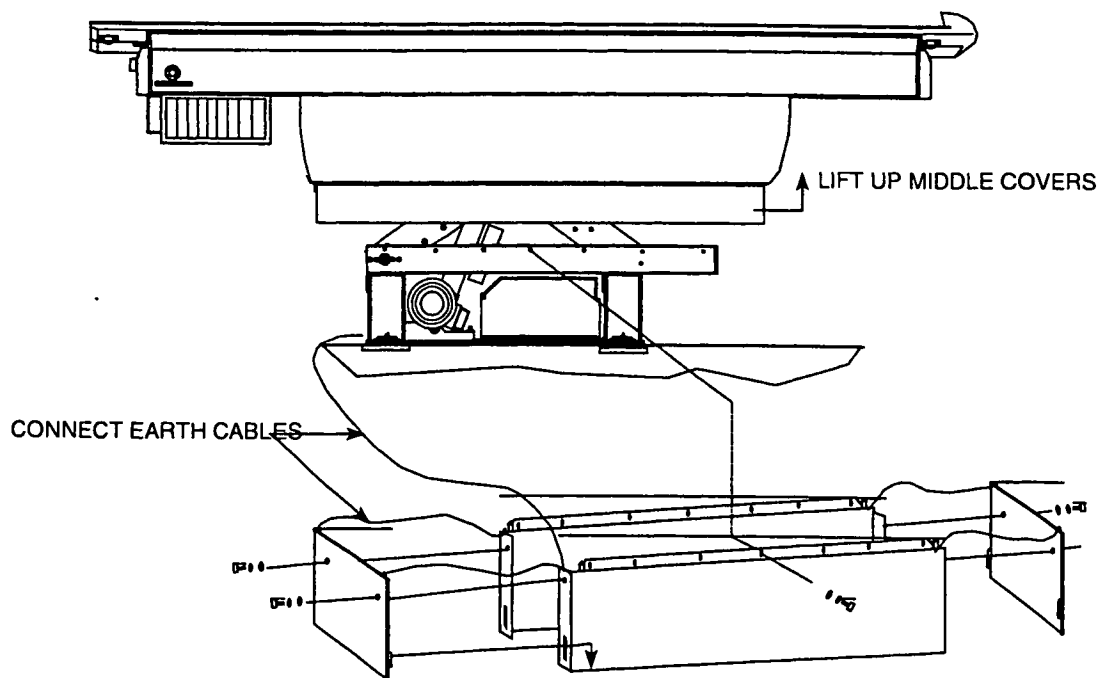

CAUTION

While moving the table to the lowest position, watch and listen if the covers are running clear of the inner frame. Adjust if necessary.

Intera 0.5 T		Intera 1.0 T		Intera 1.5 T	
Standard	Omni	Power	Omni	Power	Master

4. MOUNTING THE SYSTEM FILTER BOX COVERS

NOTE

If the system filter box back plate is mounted directly onto the finished wall, than no spacers are required. If the system filter box back plate is mounted deeper as the finished wall, than spacers are required to compensate for this extra depth


Proceed as follows:

1. Mount the top frame to the system filter box (detail C). Add spacers underneath the frame if necessary.
2. Mount the aluminium middle and upper inner covers (There is no lower inner cover).
3. Mount the upper cover. Position the cover into the top frame (detail A) and mount the bottom (detail B). Add spacers if needed.
4. Mount the middle cover. Position the middle cover over the top cover (detail B). If needed add spacers.
5. Mount the lower cover. Position the lower cover over the middle cover (detail B). Fix the lower cover with three screws on the bottom side (detail E).

Intera 0.5 T		Intera 1.0 T		Intera 1.5 T	
Standard	Omni	Power	Omni	Power	Master

Figure 24 – System Filter Box

covers

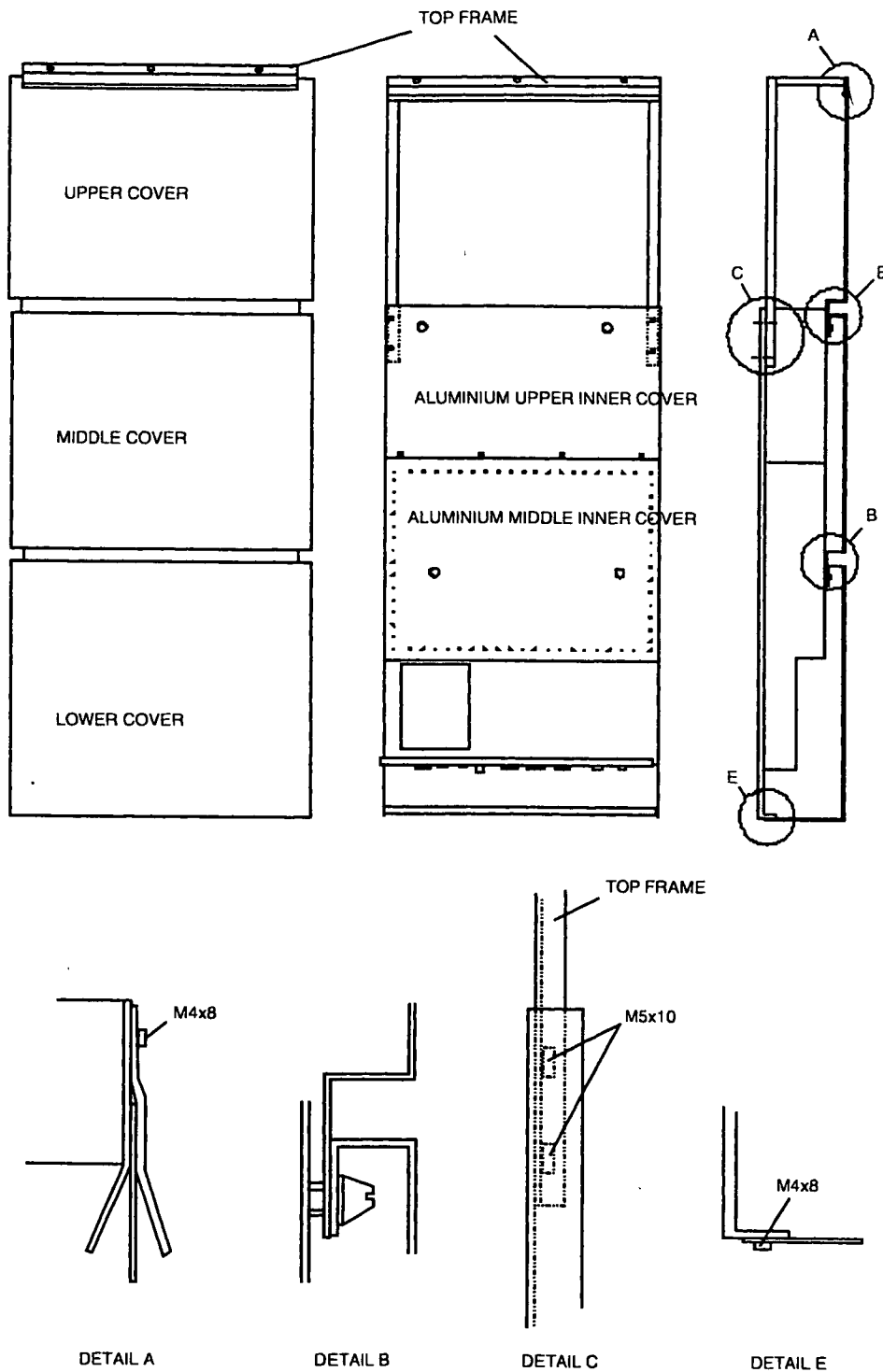

Intera 0.5 T		Intera 1.0 T		Intera 1.5 T	
Standard	Omni	Power	Omni	Power	Master

Figure 25 – System Filter Box covers

