

Intera 0 5T	Intera 1 0T		Intera 1 5T		
Standard	Omni	Power	Omni	Power	Master

Section 12

Patient Support Intera

Contents

1. SETTING TO WORK THE PATIENT SUPPORT INTERA 3

1.1. Configuration 3

1.2. Power ON of the patient support 3

2. ADJUST TABLE LENGTH MOVEMENT 4

2.1. Light-visor check 5

Intera 0.5T		Intera 1.0T		Intera 1.5T	
Standard	Omni	Power	Omni	Power	Master

Intera 0 5T		Intera 1 0T		Intera 1 5T	
Standard	Omni	Power	Omni	Power	Master

1. SETTING TO WORK THE PATIENT SUPPORT INTERA

1.1. CONFIGURATION


The configuration parameter must have been set to "MT-PATSUP" as described in section 11.

1.2. POWER ON OF THE PATIENT SUPPORT

The power-up sequence starts after you apply power to the patient support with switch F6 in the MDU. After a PLC-firmware self-check, the PLC switches into a RUN-mode. After that, the PLC firmware checks if failures are present. If a failure is detected, the TTS-led will lit for 10 seconds. After that, the PLC switches to the so-called "normal motorized mode".

Now the TTS led and the MANUAL led on the PICU will start blinking. This indicates that the table reference position must be determined. The reference switch will activate if the tabletop is about 12 mm from the mechanical end stop out. Move the tabletop by "Tumble switch" (slow movement) over this position. Both LED's will switch off.

Figure 1 – Power on self check

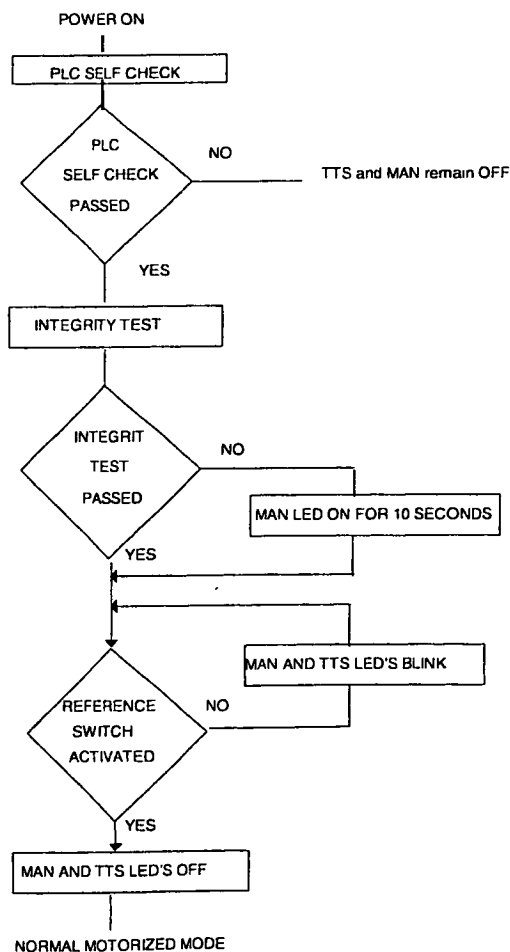

Intera 0.5T	Intera 1.0T		Intera 1.5T		
Standard	Omni	Power	Omni	Power	Master

2. ADJUST TABLE LENGTH MOVEMENT

In order to use the total length of the horizontal patient support movement, this length must be determined. Do this as follows:

- Login under TSW
- Select TPAT
- Enter: RUN ADJUSTMENTS TABLE_LENGTH

Follow the instructions on the screen. After determination of the table length, the length is stored on the host computer and into the patient support (PLC Programmable Logic Controller).

NOTE

After power OFF/ON, the value in the patient support (PLC) is lost. The 'host' computer however will compare the 'host' value against the patient support value always after start-up of the application software and at the beginning of each scan. If the values are not equal, the 'host' value will be written into the patient support (PLC).

-
- Enter: QUIT
 - Logout of TSW

Intera 0 5T		Intera 1 0T		Intera 1 5T	
Standard	Omni	Power	Omni	Power	Master

2.1. LIGHT-VISOR CHECK

The light-visor is factory mounted inside the front cone cover.

WARNING

Carefully read the safety aspects that are valid when working with a laser. Refer to section 1 chapter "Laser safety"

The distance between the light-visor and the iso-center is specified to be 827 ± 5 mm. This tolerance (± 5 mm) made it possible to build the laser into the front cone cover. This light-visor is factory adjusted.

To check the position of the laser beam, a cross is drawn at the bottom of the front cone cover. The lines have a width of 10 mm (± 5 mm). For testing, the bridge section must be removed. Check if the laser beams at the bottom of the cone-cover are within the marked cross at the bottom of the cone-cover.

Figure 2 – Laser beam position

