

Intera 0.5 T		Intera 1.0 T		Intera 1.5 T	
Standard	Omni	Power	Omni	Power	Master

Section 21.2

Service Access

Contents


1.	INTRODUCTION	3
2.	REQUIREMENTS	3
3.	INSTALLATION PROCEDURE.....	3
4.	STARTING A SERVICE SESSION IN THE GYROSCAN SOFTWARE.....	4
4.1.	Service session for a local Philips engineer.....	4
4.1.1.	Gyroscan is up and running.	4
4.1.2.	Gyroscan/TSW is not running	8
4.2.	Service session for a local NON-Philips engineer.....	9
4.2.1.	Gyroscan is up and running.	9
4.2.2.	Gyroscan/TSW is not running	10
4.3.	Service session from a remote service center.	11
5.	ENDING A SERVICE SESSION.....	12
5.1.	Gyroscan is up and running	12
5.2.	Gyroscan/TSW is not running.	13
6.	PROBLEMS IN USING THE RESPONSE GENERATOR	14

Intera 0.5 T		Intera 1.0 T		Intera 1.5 T	
Standard	Omni	Power	Omni	Power	Master

Intera 0.5 T		Intera 1.0 T		Intera 1.5 T	
Standard	Omni	Power	Omni	Power	Master

1. INTRODUCTION

A *challenge and response* mechanism, also used in the banking industry, is the basis for the Gyroscan service access control mechanism. The system will generate a challenge (string of characters) which must be responded with a string containing the correct information before protected service functionality becomes available. This response is generated by a dongle-protected program resident on the service PC. The dongle (HW key on your service PC), as already available with Philips engineers, is used to provide the authorization information which is part of the response.

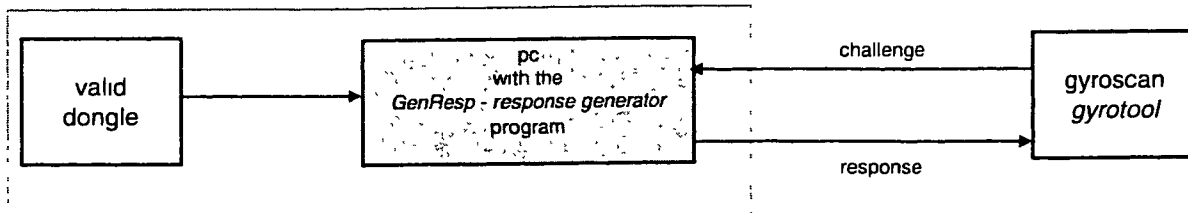

Therefor the contents of the dongle must be extended with the new MR functionality (like service level and local/remote usage).

The GenResp program, distributed with the Customer Support CD-ROM, must be installed at the engineer's service PC.

2. REQUIREMENTS

To make use of this functionality there are some requirements.

The response generator program: Version 1.3 (Check in the Help-about functions)

The response Generator program is available for the following operating systems:

OS on service PC: Windows 3.1x,
Windows 95 or Windows 98
Windows NT

The Application software on the MR systems must have one of the following releases:

Gyroscan software: R4.6.2/R5.6.2 or R6.1.2 or higher.

3. INSTALLATION PROCEDURE

Refer to the Release Bulletin for installing the response generator software on your service PC. This R.B. is delivered with the Customer Support CD-ROM and can be found in the MR section of the BBS download area 13.

Intera 0 5 T		Intera 1 0 T		Intera 1 5 T	
Standard	Omni	Power	Omni	Power	Master

4. STARTING A SERVICE SESSION IN THE GYROSCAN SOFTWARE.

Service sessions for a local Philips engineer are described in paragraph 4.1.

Service sessions for a local NON Philips engineer are described in paragraph 4.2.

Service sessions for a Remote Service Center are described in paragraph 4.3.

4.1. SERVICE SESSION FOR A LOCAL PHILIPS ENGINEER

Depending on whether or not Gyroscan is up and running, one of the following procedures can be used to obtain service access.

4.1.1 Gyroscan is up and running

4.1.2 Gyroscan/TSW is not running

4.1.1. GYROSCAN IS UP AND RUNNING.

On the console monitor place your mouse pointer in the image area and press ALT + F9 simultaneously (2 times).

Select the Manager icon.

Open a test terminal (if not already open) to get a DCL prompt (\$).

Enter: **dongle**

The next menu will appear:

```

: service functions                                11-SEP-1998 12:47
-----
1 Authorize service functions for Philips service user
2 Authorize service functions for NON-Philips service user
3 Remove authorisation(s) of service functions

-----
Select: _

```

Select: Authorize service functions for Philips service user.

The next menu will be displayed (challenge string is only an example).

```


: Authorize service funct.                          11-SEP-1998 12:57
-----
Challenge: 56xdp3k5
Response  _____
User name  _____

-----
1-----1
«Cancel» or «Proceed» to continue


```

Intera 0.5 T		Intera 1.0 T		Intera 1.5 T		
Standard	Omni	Power	Omni	Power	Master	


On your service PC, connect the dongle HW key to the parallel port.
 Startup Windows and select the Program Group "PMS Service Workbench"
 Double click on the Response generator icon.

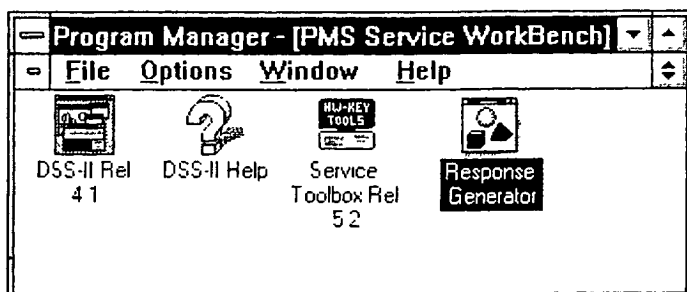

A hardware key password dialog will be opened.

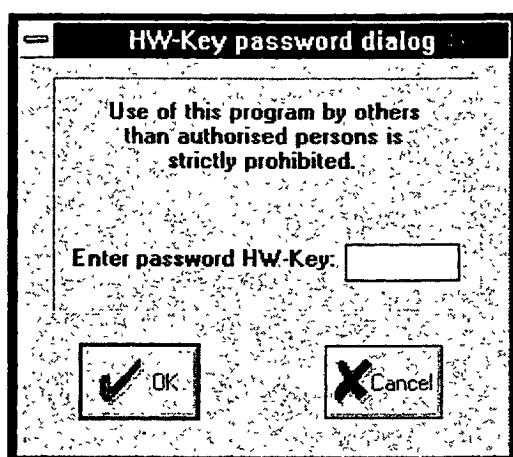

Fill in the password of your HW key.
 Select: <OK>
 The response generator will be opened.

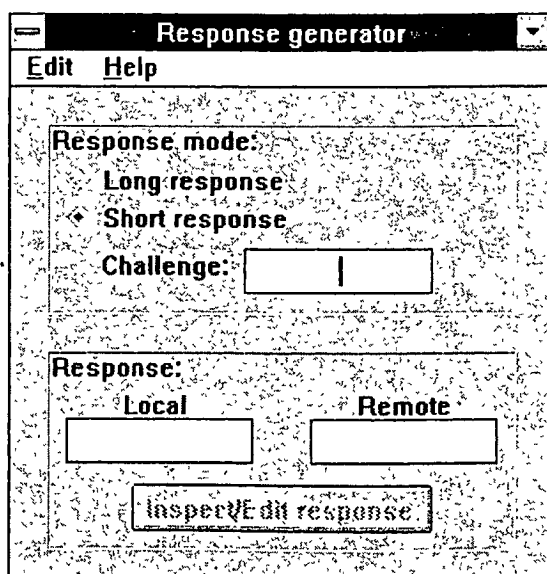

Intera 0.5 T		Intera 1.0 T		Intera 1.5 T	
Standard	Omni	Power	Omni	Power	Master

Enter the Challenge string on your response generator.

Response generator

Edit Help

Response mode:

- > Long response
- ◆ Short response

Challenge: 56xdp3k5

Response:

Local Remote

75v5-6zqy 4td9-e6c6

Inspect/Edit response

A response will be generated.

If required you can change some items like the authorization level, with the <Inspect/Edit response> button and generate a different response with the <OK> button.

Short response information

Dongle identification:

Access mode

☒ Local ☐ Remote

Enabled authorisation levels

☒ 0 ☒ 1 ☐ 2 ☐ 3

OK Cancel

Enter the local response in your Gyroscan SW.

Fill in your name and press the <Proceed> key on your keyboard.

Intera 0.5 T		Intera 1.0 T		Intera 1.5 T	
Standard	Omni	Power	Omni	Power	Master

· Authorize service funct	11-SEP-1998 13:15

Challenge: 56xdp3k5	
Response	75v5-6zqy
User name	j. janssen_____
1-----1	
<Cancel> or <Proceed> to continue	

A service software key will now be generated and will be placed on position 23 or 24 of the temporary keys.

To quit the dongle application press the <Cancel> key on your keyboard.

Go back to your Gyroscan window by double clicking on the Gyroview icon.

The system is ready for using the service functions for which you are authorized by your HW key.

The first time a "service" window is entered the next message will be displayed:

You enter a service protected function.

Performing checks and/or adjustments and/or modifying any of the involved parameters affects the imaging capabilities of the Gyroscan system.
This may only be performed by persons licensed by Philips Medical Systems.

Press <Proceed> to enter.

Intera 05 T		Intera 10 T		Intera 15 T	
Standard	Omni	Power	Omni	Power	Master

4.1.2. GYROSCAN/TSW IS NOT RUNNING

In case you want to use protected ASW or TSW functionality you must authorize first the service functions using Gyrotool. In case the Gyroscan application is not up and running, you can start your service session also with:

Login : **Gyrotool**

Select: Modify system configuration

Select: Options (SW Keys)

```

Figman      : Options (SW-keys)          GYRO$SITE:gyrofig.idx; 11-SEP-1998 13:35
-----
 1 Display / Register permanent key
 2 View SW options, modify option status
 3 Display / Register temporary keys
 4 Display options authorized by keys
 5 Display keys which authorize options
 6 Service options
 7 Update locktype

1-----1
Select: _

```

Select: Service options

```

Figman      : service functions          GYRO$SITE:gyrofig.idx; 11-SEP-1998 13:36
-----
 1 Authorize service functions for Philips service user
 2 Authorize service functions for NON-Philips service user
 3 Remove authorisation(s) of service functions

-----
Select: _

```

Select: Authorize service functions for Philips service user.

The rest of the procedure is the same as described in paragraph 4.1.1

Intera 0.5 T		Intera 1.0 T		Intera 1.5 T	
Standard	Omni	Power	Omni	Power	Master

4.2. SERVICE SESSION FOR A LOCAL NON-PHILIPS ENGINEER.

For a non-Philips engineer a so called long response is used which authorizes him or her to do service work on 1 specific system for a specified time. This long response can be obtained at the Customer Support department of Philips Medical Systems (Refer to Newsletter MR-105).

Also here on of the following procedures can be used:

- 4.2.1 Gyroscan is up and running
- 4.2.2 Gyroscan/TSW is not running

4.2.1. GYROSCAN IS UP AND RUNNING.

On the console monitor place your mouse pointer in the image area and press ALT + F9 simultaneously.

Double click on the Manager icon.

Open a test terminal to get a DCL prompt.

Enter: **dongle**

```

: service functions                                     11-SEP-1998 12:47
-----
1 Authorize service functions for Philips service user
2 Authorize service functions for NON-Philips service user
3 Remove authorisation(s) of service functions
-----
Select: _

```

Select: Authorize service functions for NON-Philips service user.

The next menu will be displayed.

Fill in the long response and your name and press <Proceed>

```

Figman      : Authorize service funct.  GYRO$SITE:gyrofig.idx; 11-SEP-1998 14:12
-----
Response    de9u-5gne-q2la-eyeb
User name    f.hals_____
-----
1-----1
«Cancel» or «Proceed» to continue

```

A service software key will now be generated and will be placed on position 23 or 24 of the temporary keys.

To quit the dongle application select: <Cancel>

Go back to your Gyroscan window by double clicking on the Gyroview icon.

The system is ready for using the service functions for which you are authorized

The first time a 'service' window is entered the next message will be displayed:

```

You enter a service protected function.

Performing checks and/or adjustments and/or modifying any of the involved
parameters affects the imaging capabilities of the gyroscan system.
This may only be performed by persons licensed by Philips Medical Systems.

```

Press <Proceed> to enter.

Intera 0.5 T		Intera 1.0 T		Intera 1.5 T	
Standard	Omni	Power	Omni	Power	Master

4.2.2. GYROSCAN/TSW IS NOT RUNNING

When the Gyroscan or TSW application is not up and running, you can start your service session with:

Login **Gyrottool**

Select: Modify system configuration

Select: Options (SW Keys).

```

Figman      : Options (SW-keys)          GYROSSITE:gyrofig.idx; 11-SEP-1998 13:35
-----
1  Display / Register permanent key
2  View SW options, modify option status
3  Display / Register temporary keys
4  Display options authorized by keys
5  Display keys which authorize options
6  Service options
7  Update locktype

1-----1
Select: _

```

Select: Service options

```

Figman      : service functions          GYROSSITE:gyrofig.idx; 11-SEP-1998 13:36
-----
1  Authorize service functions for Philips service user
2  Authorize service functions for NON-Philips service user
3  Remove authorisation(s) of service functions

-----
Select: _

```

The rest of the procedure is the same as described in paragraph 4.2.1

Intera 0.5 T		Intera 1.0 T		Intera 1.5 T	
Standard	Omni	Power	Omni	Power	Master

4.3. SERVICE SESSION FROM A REMOTE SERVICE CENTER.

You login remotely to the system via the User : **MODEM**

Enter: **\$gyroscan**

Enter: Gyrotest>**dongle**

The rest of the procedure is the same as described in paragraph 4.1.1 except that you can use the cut and paste mechanism for entering the strings.

Response generator

Edit Help

Response mode:

☐ Long response

☒ Short response

Challenge: 56xdp3k5

Response:

Local: 75v5-6zqy Remote: 4td9-e6c6

Inspect/Edit response

The Response Generator program has an edit mode where it is possible to copy (CTRL + C) the string from the remote output block and paste (CTRL + V) it to the response string in the dongle program.

: Authorize service funct. 11-SEP-1998 13:15

Challenge: 56xdp3k5

Response 4td9-e6c6

User name RTAC

1-----1

<Cancel> or <Proceed> to continue

Intera 0 5 T		Intera 1 0 T		Intera 1 5 T	
Standard	Omni	Power	Omni	Power	Master

5. ENDING A SERVICE SESSION.

At the end of a service session the service user needs to remove the service key, to prevent that other unauthorized persons make use of the service functions.

The key generated with the short response will be valid till the end of the next day: at 24:00 hours it will expire automatically.

The key generated with the long response will normally be valid for a longer time.

Therefore we advise to remove the keys when finishing the service session.

There are two ways to remove the service keys. The procedure is the same for the Philips user and the NON-Philips users and can be done also from a Remote Service Center.

5.1. GYROSCAN IS UP AND RUNNING

Put your mouse pointer in the image area and press ALT + F9 simultaneously.

Double click on the Manager icon (if needed).

Open a test terminal (if not already opened) to get a DCL prompt (\$).

Enter: **dongle**

```

: service functions                                     11-SEP-1998 12:47
-----
1  Authorize service functions for Philips service user
2  Authorize service functions for NON-Philips service user
3  Remove authorisation(s) of service functions

-----
Select: _

```

Select: Remove authorization(s) of service functions.

```

Flgman      :remove service functions  GYRO$SITE:gyrofig.idx; 11-SEP-1998 15:28
-----
identification

1 10100 15-SEP-1999 loc oth 0,1 f.hals

-----
Select :

```

Select the authorization which you want to remove (not the key will be displayed but the user information).

Press: **<Return>**

The service key will be removed.

Press the **<Cancel>** key on the keyboard to leave the dongle program.

At certain points in the ASW/TSW functions the SW keys are checked. If you pass one of these points and try to enter a protected function when the key is removed you will not be able to enter these functions.

Intera 0.5 T		Intera 1.0 T		Intera 1.5 T	
Standard	Omni	Power	Omni	Power	Master

5.2. GYROSCAN/TSW IS NOT RUNNING.

After a TSW session or when the Gyroscan application is not up and running, you can end your service session with:

Login · **Gyrottool**

Select: Modify system configuration

Select: Options (SW Keys)

```

Figman      : Options (SW-keys)      GYRO$SITE:gyrofig.idx; 11-SEP-1998 13:35
-----
1 Display / Register permanent key
2 View SW options, modify option status
3 Display / Register temporary keys
4 Display options authorized by keys
5 Display keys which authorize options
6 Service options
7 Update locktype

1-----1
Select: _

```

Select: Service options

```

Figman      : service functions      GYRO$SITE:gyrofig.idx; 11-SEP-1998 13:36
-----
1 Authorize service functions for Philips service user
2 Authorize service functions for NON-Philips service user
3 Remove authorisation(s) of service functions

-----
Select: _

```

Select: Remove authorization(s) of service functions.

```

Figman      :remove service functions  GYRO$SITE:gyrofig.idx; 11-SEP-1998 15:28
-----
identification
1 10100 15-SEP-1999 loc oth 0,1 f.hals

-----
Select :

```

Select the authorization to remove.

Press: <Return>

The service key will be removed.

Press the <Cancel> key on the keyboard to leave the dongle program.

Intera 0.5 T		Intera 1.0 T		Intera 1.5 T	
Standard	Omni	Power	Omni	Power	Master

6. PROBLEMS IN USING THE RESPONSE GENERATOR

- **Error during startup of the response generator program.**

When you get the following error after starting the program:

'HW key was not detected. Check for proper connection'

Solution : Check that you have installed the HW key to the parallel port of the service PC. Close your window application and startup windows again.

- **Error during startup of the response generator program**

When you get the following error after starting the program:

'Your key is not authorized for usage of this software'

Solution : You have a key connected which has not been update to use the response generator program. Ask for a updated key contents via the BBS (Refer to Newsletter MR-105). In the mean time you can ask the RS Center to generate a valid string for the system.